


CENGAGE LEARNING LABORATORY SERIES *for*
Organic Chemistry

A Microscale Approach *to* Organic Laboratory Techniques

SIXTH EDITION

Donald L. Pavia
Gary M. Lampman
George S. Kriz
Randall G. Engel

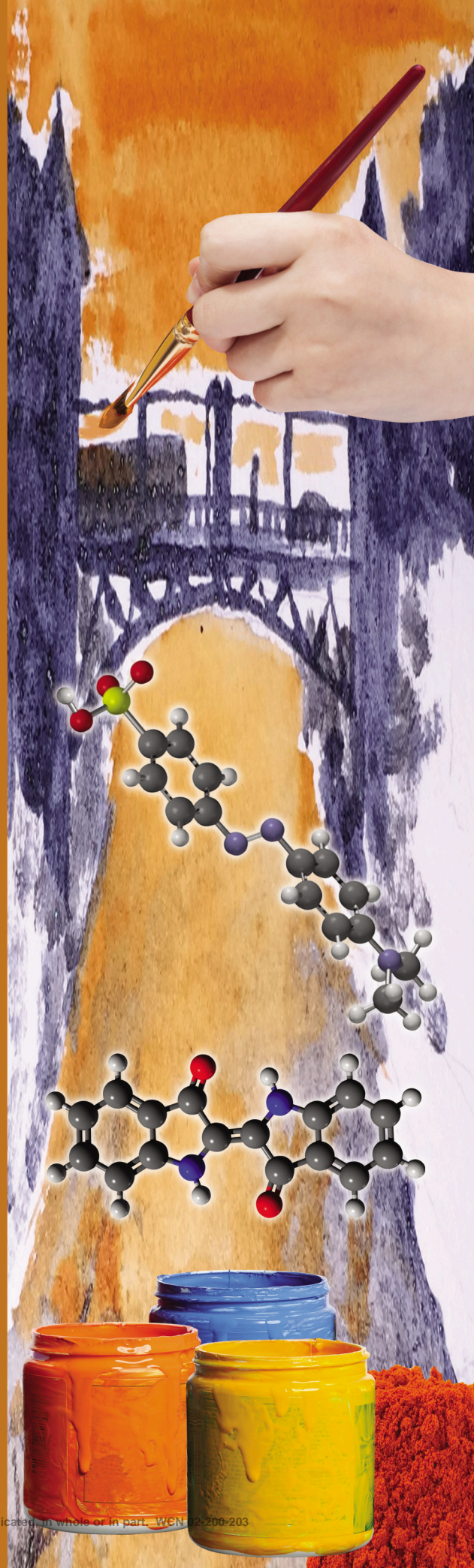

A Microscale
Approach
to
Organic
Laboratory
Techniques

SIXTH EDITION

A Microscale Approach *to* Organic Laboratory Techniques

SIXTH EDITION

Donald L. Pavia
Gary M. Lampman
George S. Kriz

Western Washington University
Bellingham, Washington

Randall G. Engel

North Seattle Community College
Seattle, Washington



Australia • Brazil • Japan • Korea • Mexico • Singapore • Spain • United Kingdom • United States

***A Microscale Approach to Organic
Laboratory Techniques, Sixth Edition***
Donald L. Pavia, George S. Kriz, Gary M.
Lampman, and Randall G. Engel

Product Director: Dawn Giovannello

Associate Product Manager: Courtney Heilman

Content Developer: Brendan R Killion

Product Assistant: Kristina Cannon

Marketing Manager: Janet del Mundo

Art and Cover Direction, Production
Management, and Composition: Lumina
Datamatics, Inc.

Manufacturing Planner: Judy Inouye

Cover Image: R. Gino Santa maria/Shutterfree,
Lic./Dreamstime.com; © Petr Vodicka |
Dreamstime.com; vvoe/Fotolia LLC;
marylooo/iStockphoto; © Donald Pavia;
© Ailish O'Sullivan

Unless otherwise noted all items

© Cengage Learning®

© 2018, 2013 Cengage Learning

ALL RIGHTS RESERVED. No part of this work covered by the copyright herein may be reproduced or distributed in any form or by any means, except as permitted by U.S. copyright law, without the prior written permission of the copyright owner.

For product information and technology assistance, contact us at
Cengage Learning Customer & Sales Support, 1-800-354-9706.

For permission to use material from this text or product,
submit all requests online at **www.cengage.com/permissions.**

Further permissions questions can be e-mailed to
permissionrequest@cengage.com.

Library of Congress Control Number: 2016951799

Student Edition:

ISBN: 978-1-305-96834-9

Cengage Learning

20 Channel Center Street
Boston, MA 02210
USA

Cengage Learning is a leading provider of customized learning solutions with employees residing in nearly 40 different countries and sales in more than 125 countries around the world. Find your local representative at **www.cengage.com.**

Cengage Learning products are represented in Canada by
Nelson Education, Ltd.

To learn more about Cengage Learning Solutions, visit **www.cengage.com.**

Purchase any of our products at your local college store or at our preferred online store **www.cengagebrain.com.**

*This book is dedicated to
our organic chemistry laboratory students*

Preface

STATEMENT OF MISSION AND PURPOSE IN REVISING THE TEXTBOOK

The purpose of this lab book is to teach students the techniques of organic chemistry. We desire to share our love of the organic chemistry lab and the joy it brings us with our students! In this edition, we have provided many new, up-to-date experiments that will demonstrate how organic chemistry is evolving. We have updated and improved many of the standard experiments from previous editions, and we have added some new ones. For example, we have included some experiments involving dyes and soap. To make the connection of organic chemistry to our everyday world even more real, we have added a project experiment that asks the students to formulate a paint and then use it in an art project. We think that you will be enthusiastic about this new edition. Many of the new experiments will not be found in other laboratory manuals, but we have been careful to retain all of the standard reactions and techniques, such as the Friedel-Crafts reaction, aldol condensation, Grignard synthesis, and basic experiments designed to teach crystallization, chromatography, and distillation.

SCALE IN THE ORGANIC LABORATORY

Experiments in organic chemistry can be conducted at different scales using varying amounts of chemicals and different styles of glassware. We have two versions of our laboratory textbooks that teach organic laboratory techniques. Our microscale book (*A Microscale Approach to Organic Laboratory Techniques, Sixth Edition*) makes use of 14/10 standard tapered glassware. Our version of a “macroscale” textbook (*A Small Scale Approach to Organic Laboratory Techniques, Fourth Edition*) uses the traditional larger scale 19/22 standard tapered glassware. The fourth edition of our small scale book was published in 2016.

Over the years that we have been involved with developing experiments, we have learned that students can easily adjust to working with the small laboratory equipment that is used in this microscale book. As students and faculty learn to appreciate the impact of laboratory classroom experiments on the environment, they become more aware that it is not necessary to consume large quantities of chemicals. Students come to appreciate the importance of reducing waste generated in the organic laboratory. All of us, students and faculty alike, are becoming more “green.”

MAJOR FEATURES OF THE TEXTBOOK THAT WILL BENEFIT THE STUDENT

When we published our first organic laboratory textbook in 1976, a major goal was to demonstrate to students how organic chemistry significantly impacts our lives in the real world. This was accomplished by including experiments with a real-world connection and by including many topical essays that related the experiments to everyday world applications. In this edition, we have taken this emphasis to a new level. For example, we have added two new experiments involving the synthesis of two widely used dyes, methyl orange and indigo. These dyes can then be used to formulate a paint in the experiment Formulation of a Paint and Art Project. Not only do students learn about the chemistry involved in the formulation of a paint, but they also paint a picture of their own creation. Many students at North Seattle College and the University of Washington report that this is one of their favorite experiments in the organic laboratory class! We have also added a new essay on Dyes that gives further examples of how these new experiments are related to our everyday lives.

Another real-world experiment that we are especially excited about is Preparation of Soap. This experiment was developed by one of our organic chemistry students, who is a professional soap maker! Students learn about the chemistry of soap making, and they make a bar of soap that can be used at home. We have also included a new essay on Soap.

A number of experiments are linked together to create multistep syntheses. The advantage of this approach is that you will be doing something different from your neighbor in the laboratory. Wouldn't you like to be carrying out an experiment that is not the same as your neighbor's? Maybe you will be synthesizing a new compound that hasn't been reported in the chemical literature! You and your fellow students will not all be doing the same reaction on the same compound: for example, some of you will be carrying out the chalcone reaction, others the "green" epoxidation, and still others the cyclopropanation of the resulting chalcones.

GREEN CHEMISTRY

We have continued an emphasis on Green Chemistry in this edition. The Green Chemistry experiments decrease the need for hazardous waste disposal, leading to reduced contamination of the environment. These experiments use less toxic reactants and solvents. For example, water is used as a solvent in some experiments. Almost all experiments have been reduced in scale compared to the traditional macroscale experiments. Experiments that are particularly good for illustrating the Green Chemistry approach include Biodiesel, Chiral Reduction of Ethyl Acetoacetate, Aqueous-Based Organozinc Reactions, Grubbs-Catalyzed Metathesis of Eugenol with 1,4-Butanediol, Diels-Alder Reaction with Anthracene-9-methanol, and Green Epoxidation of Chalcones. We have also added a new Green oxidation reaction using Oxone® in an Oxidation-Reduction Scheme: Borneol, Camphor, Isoborneol. Oxone® is a more reliable alternative to bleach, which we have used in previous editions of this textbook.

In keeping with the Green Chemistry approach, we have suggested an alternative way of approaching qualitative analysis. This approach makes extensive use of spectroscopy to solve the structure of organic unknowns. In this approach, some of the traditional tests have been retained, but the main emphasis is on using

spectroscopy. In this way, we have attempted to show students how to solve structures in a more modern way, similar to that used in a research laboratory. The added advantage to this approach is that waste is considerably reduced.

NEW TO THIS EDITION

Many of the new experiments in this edition demonstrate the relationship between organic chemistry and our everyday lives. This edition also includes updating of the essays and the chapters on techniques. New experiments added for this edition include:

- Experiment 26 Preparation of Soap
- Experiment 33 An Oxidation-Reduction Scheme: Borneol, Camphor, Isoborneol
- Experiment 46 Preparation of Methyl Orange
- Experiment 47 Preparation of Indigo
- Experiment 48 Formulation of a Paint and Art Project

New Essays include:

- Soap
- Dyes

As in previous editions, the technique chapters include both microscale and macroscale techniques. Many of the references in the technique chapters have been updated. New material on diastereotopic protons has been added to Technique 26, Nuclear Magnetic Resonance Spectroscopy. Technique 29, Guide to the Chemical Literature, has been revised.

SUPPORTING RESOURCES

Please visit <http://www.cengage.com/chemistry/pavia/microorglab6e> for information about student and instructor resources for this text.

ACKNOWLEDGMENTS

We owe our sincere thanks to the many colleagues who have used our textbooks and who have offered their suggestions for changes and improvements to our laboratory procedures or discussions. Although we cannot mention everyone who has made important contributions, we must make special mention of Albert Burns (North Seattle College), Charles Wandler (Western Washington University), Emily Borda (Western Washington University), Frank Deering (North Seattle College), Jacob Frank (North Seattle College), Gregory O'Neil (Western Washington University), James Vyvyan (Western Washington University), Khushroo Daruwala (University of Washington Bothell), Scott Clary (North Seattle College), and Timothy Clark (University of San Diego).

In preparing this new edition, we have also attempted to incorporate the many improvements and suggestions that have been forwarded to us by the many instructors who have been using our materials over the past several years.

We are especially grateful to James Patterson, faculty member of North Seattle College, who has given us permission to include several of his experiments in our

textbooks. His ideas and enthusiastic support of our textbooks for many years have contributed immensely to the success of our textbooks.

We thank all who contributed, with special thanks to our Senior Product Manager, Lisa Lockwood; Associate Content Developer, Brendan Killion; Content Project Manager, James Zayicek; Associate Marketing Manager, Ana Albinson; and Associate Program Manager, Sharib Asrar at Lumina Datamatics.

We are especially grateful to the students and friends who have volunteered to participate in the development of experiments or who offered their help and criticism. We owe special thanks to Sean Ichiun Choe, organic chemistry student at North Seattle College, who developed and wrote most of Experiment 24 (Preparation of Soap). Sean's expertise as a soap maker in the real world is reflected in this valuable addition to our book. Sean also made valuable contributions to the Soap essay.

We are also grateful to Alish O'Sullivan, student at North Seattle College, who painted the picture of the Montlake Bridge, which appears on the cover of this textbook. This painting was created by Alish while performing the new experiment, Formulation of a Paint and Art Project, which appears in this textbook.

Finally, we wish to thank our families and special friends, especially Neva-Jean Pavia, Marian Lampman, and Karin Granstrom, for their encouragement, support, and patience.

Donald L. Pavia (pavia@chem.wvu.edu)
Gary M. Lampman (lampman@chem.wvu.edu)
George S. Kriz (George.Kriz@wvu.edu)
Randall G. Engel (Randall.Engel@seattlecolleges.edu)

August 2016

How To Use This Book

OVERALL STRUCTURE OF THE BOOK

This textbook is divided into two major sections (see Table of Contents). The first section, which includes Part One through Part Five, contains all of the experiments in this book. The second major section includes only Part Six, which contains all of the important techniques you will use in performing the experiments in this book. Interspersed among the experiments in Part One through Part Three is a series of essays. The essays provide a context for many of the experiments and often relate the experiment to real world applications. When your instructor assigns an experiment, he or she will often assign an essay and/or several techniques chapters along with the experiment. Before you come to lab, you should read all of these. In addition, it is likely that you will need to prepare some sections in your laboratory notebook (see Technique 2) before you come to the lab.

STRUCTURE OF THE EXPERIMENTS

In this section we discuss how each experiment is organized in the textbook. To follow this discussion, you may want to refer to a specific experiment, such as Experiment 13.

Multiple Parts Experiments

Some experiments, such as Experiment 13, are divided into two or more individual parts that are designated by the experiment number and the letters A, B, etc. In some experiments, like Experiment 13, each part is a separate but related experiment, and you will most likely perform only one part. In Experiment 13, you would do Experiment 13A (Isolation of Caffeine from Tea Leaves) or Experiment 13B (Isolation of Caffeine from a Tea Bag). In other experiments, for example Experiment 32, the various parts can be linked together to form a multistep synthesis. In a few experiments, such as Experiment 22, the last part describes how you should analyze your final product.

Featured Topics and Techniques Lists

Directly under the title of each experiment (see Experiment 13), there will be a list of topics. These topics may explain what kind of experiment it is, such as isolation of a natural product or Green Chemistry. The topics may also include major techniques that are required to perform the experiment, such as crystallization or extraction.

Required Reading

In the introduction to each experiment, there will be a section labeled Required Reading. Within this section, some of the required readings are labeled Review and some are labeled New. You should always read the chapters listed in the New section. Sometimes it will also be helpful to do the readings in the Review section.

Special Instructions

You should always read this section since it may include instructions that are essential to the success of the experiment.

Suggested Waste Disposal

This very important section gives instructions on how to dispose of the waste generated in an experiment. Often your instructor will provide you with additional instructions on how to handle the waste.

Notes to Instructor

It will usually not be necessary to read this section. This section provides special advice for the instructor that will help to make the experiment successful.

Procedure

This section provides detailed instructions on how to carry out the experiments. Within the procedure, there will be many references to the techniques chapters, which you may need to consult in order to perform an experiment.

Report

In some experiments, specific suggestions for what should be included in the laboratory report will be given. Your instructor may refer to these recommendations or may have other directions for you to follow.

Questions

At the end of most experiments will be a list of questions related to the experiment. It is likely that your instructor will assign at least some of these questions along with the laboratory report.

Contents

Preface vii

PART 1 Introduction to Basic Laboratory Techniques 1

- 1 Introduction to Microscale Laboratory 2
- 2 Solubility 12
- 3 Crystallization 20
 - 3A Semimicroscale Crystallization—Erlenmeyer Flask and Hirsch Funnel 21
 - 3B Microscale Crystallization—Craig Tube 24
 - 3C Selecting a Solvent to Crystallize a Substance 26
 - 3D Mixture Melting Points 27
 - 3E Critical Thinking Application 28
- 4 Extraction 32
 - 4A Extraction of Caffeine 33
 - 4B Distribution of a Solute between Two Immiscible Solvents 35
 - 4C How Do You Determine Which One Is the Organic Layer? 36
 - 4D Use of Extraction to Isolate a Neutral Compound from a Mixture Containing an Acid or Base Impurity 37
 - 4E Critical Thinking Application 39
- 5 A Separation and Purification Scheme 42
- 6 Chromatography 45
 - 6A Thin-Layer Chromatography 46
 - 6B Selecting the Correct Solvent for Thin-Layer Chromatography 48
 - 6C Monitoring a Reaction with Thin-Layer Chromatography 49
 - 6D Column Chromatography 50
- 7 Infrared Spectroscopy and Boiling-Point Determination 54
- 8 Simple and Fractional Distillation 58
 - 8A Simple and Fractional Distillation (Semimicroscale Procedure) 60
 - 8B Simple and Fractional Distillation (Microscale Procedure) 64
- Essay Aspirin 66
 - 9 Acetylsalicylic Acid 69
- Essay Analgesics 73
- 10 Isolation of the Active Ingredient in an Analgesic Drug 77
- 11 Acetaminophen 81
 - 11A Acetaminophen (Microscale Procedure) 82
 - 11B Acetaminophen (Semimicroscale Procedure) 84
- Essay Identification of Drugs 87
- 12 TLC Analysis of Analgesic Drugs 89

Essay	Caffeine	94
13	Isolation of Caffeine from Tea or Coffee	98
	13A Extraction of Caffeine from Tea with Methylene Chloride	101
	13B Extraction of Caffeine from Tea or Coffee Using Solid Phase Extraction (SPE)	103
Essay	Esters—Flavors and Fragrances	107
14	Isopentyl Acetate (Banana Oil)	110
	14A Isopentyl Acetate (Microscale Procedure)	111
	14B Isopentyl Acetate (Semimicroscale Procedure)	113
Essay	Terpenes and Phenylpropanoids	116
15	Essential Oils: Extraction of Oil of Cloves by Steam Distillation	120
	15A Oil of Cloves (Microscale Procedure)	121
	15B Oil of Cloves (Semimicroscale Procedure)	123
Essay	Stereochemical Theory of Odor	125
16	Spearmint and Caraway Oil: (+)- and (–)-Carvones	129
Essay	The Chemistry of Vision	137
17	Isolation of Chlorophyll and Carotenoid Pigments from Spinach	142
Essay	Ethanol and Fermentation Chemistry	149
18	Ethanol from Sucrose	152
PART 2	Introduction to Molecular Modeling	157
Essay	Molecular Modeling and Molecular Mechanics	158
19	An Introduction to Molecular Modeling	163
	19A The Conformations of <i>n</i> -Butane: Local Minima	164
	19B Cyclohexane Chair and Boat Conformations	165
	19C Substituted Cyclohexane Rings (Critical Thinking Exercises)	166
	19D <i>cis</i> - and <i>trans</i> -2-Butene	166
Essay	Computational Chemistry— <i>ab Initio</i> and Semiempirical Methods	168
20	Computational Chemistry	176
	20A Heats of Formation: Isomerism, Tautomerism, and Regioselectivity	177
	20B Heats of Reaction: S _N 1 Reaction Rates	178
	20C Density–Electrostatic Potential Maps: Acidities of Carboxylic Acids	179
	20D Density–Electrostatic Potential Maps: Carbocations	180
	20E Density–LUMO Maps: Reactivities of Carbonyl Groups	180
PART 3	Properties and Reactions of Organic Compounds	183
21	Reactivities of Alkyl Halides	184
22	Nucleophilic Substitution Reactions: Competing Nucleophiles	189
	22A Competitive Nucleophiles with 1-Butanol or 2-Butanol	191
	22B Competitive Nucleophiles with 2-Methyl-2-Propanol	193
	22C Analysis	194
23	Synthesis of <i>n</i> -Butyl Bromide and <i>t</i> -Pentyl Chloride	198
	23A <i>n</i> -Butyl Bromide	200
	23B <i>n</i> -Butyl Bromide (Semimicroscale Procedure)	202
	23C <i>t</i> -Pentyl Chloride (Microscale Procedure)	203
	23D <i>t</i> -Pentyl Chloride (Semimicroscale Procedure)	204
	23E <i>t</i> -Pentyl Chloride (Macroscale Procedure)	205
24	4-Methylcyclohexene	207
	24A 4-Methylcyclohexene (Microscale Procedure)	209
	24B 4-Methylcyclohexene (Semimicroscale Procedure)	210

Essay	Fats and Oils	213
25	Methyl Stearate from Methyl Oleate	218
Essay	Soap	223
26	Preparation of Soap	227
	26A Preparation of Soap from 70% Lard and 30% Coconut Oil	229
	26B Preparation of Several Soaps with a Given % Composition	231
Essay	Petroleum and Fossil Fuels	234
27	Gas-Chromatographic Analysis of Gasolines	243
Essay	Biofuels	248
28	Biodiesel	252
	28A Biodiesel from Coconut Oil	254
	28B Biodiesel from Other Oils	255
	28C Analysis of Biodiesel	255
Essay	Green Chemistry	258
29	Chiral Reduction of Ethyl Acetoacetate; Optical Purity Determination	264
	29A Chiral Reduction of Ethyl Acetoacetate	265
	29B NMR Determination of the Optical Purity of Ethyl (S)-3-Hydroxybutanoate	269
30	Nitration of Aromatic Compounds Using a Recyclable Catalyst	274
31	Reduction of Ketones Using Carrots as Biological Reducing Agents	278
32	Resolution of (\pm)- α -Phenylethylamine and Determination of Optical Purity	281
	32A Resolution of (\pm)- α -Phenylethylamine	283
	32B Determination of Optical Purity Using NMR and a Chiral Resolving Agent	287
33	An Oxidation–Reduction Scheme: Borneol, Camphor, Isoborneol	289
34	Multistep Reaction Sequences: The Conversion of Benzaldehyde to Benzilic Acid	304
	34A Preparation of Benzoin by Thiamine Catalysis	305
	34B Preparation of Benzil	311
	34C Preparation of Benzilic Acid	313
35	Triphenylmethanol and Benzoic Acid	317
	35A Triphenylmethanol	322
	35B Benzoic Acid	324
36	Aqueous-Based Organozinc Reactions	328
37	Sonogashira Coupling of Iodosubstituted Aromatic Compounds with Alkynes using a Palladium Catalyst	332
38	Grubbs-Catalyzed Metathesis of Eugenol with 1,4-Butenediol to Prepare a Natural Product	342
39	The Aldol Condensation Reaction: Preparation of Benzalacetophenones (Chalcones)	349
40	A Green Enantioselective Aldol Condensation Reaction	354
41	Preparation of an α,β -Unsaturated Ketone via Michael and Aldol Condensation Reactions	361
42	Preparation of Triphenylpyridine	366
43	The Wittig Reaction: Preparation of 1,4-Diphenyl-1,3-butadiene	369
	43A Benzyltriphenylphosphonium Chloride (Wittig Salt)	372
	43B Preparation of 1,4-Diphenyl-1,3-Butadiene Using Sodium Ethoxide to Generate the Ylide	372
	43C Preparation of 1,4-Diphenyl-1,3-Butadiene Using Potassium Phosphate to Generate the Ylide	374

44	Relative Reactivities of Several Aromatic Compounds	377
45	Nitration of Methyl Benzoate	381
	Essay Synthetic Dyes	386
46	Preparation of Methyl Orange	392
47	Preparation of Indigo	396
48	Formulation of a Paint and Art Project	399
	Essay Local Anesthetics	402
49	Benzocaine	406
	Essay Pheromones: Insect Attractants and Repellents	410
50	<i>N,N</i> -Diethyl- <i>m</i> -toluamide: The Insect Repellent “OFF”	418
	Essay Sulfa Drugs	423
51	Sulfa Drugs: Preparation of Sulfanilamide	426
	Essay Polymers and Plastics	431
52	Preparation and Properties of Polymers: Polyester, Nylon, and Polystyrene	441
	52A Polyesters	442
	52B Polyamide (Nylon)	443
	52C Polystyrene	445
	52D Infrared Spectra of Polymer Samples	446
	Essay Diels–Alder Reaction and Insecticides	449
53	The Diels–Alder Reaction of Cyclopentadiene with Maleic Anhydride	455
54	The Diels–Alder Reaction with Anthracene-9-methanol	459
55	Photoreduction of Benzophenone and Rearrangement of Benzopinacol to Benzopinacolone	462
	55A Photoreduction of Benzophenone	463
	55B Synthesis of β -Benzopinacolone: The Acid-Catalyzed Rearrangement of Benzopinacol	469
	Essay Fireflies and Photochemistry	471
56	Luminol	474
PART 4	Identification of Organic Substances	479
57	Identification of Unknowns	480
	57A Solubility Tests	487
	57B Tests for the Elements (N, S, X)	494
	57C Tests for Unsaturation	500
	57D Aldehydes and Ketones	504
	57E Carboxylic Acids	510
	57F Phenols	512
	57G Amines	515
	57H Alcohols	519
	57I Esters	523
PART 5	Project-Based Experiments	501
58	Preparation of a C-4 or C-5 Acetate Ester	528
59	Competing Nucleophiles in S_N1 and S_N2 Reactions: Investigations Using 2-Pentanol and 3-Pentanol	532
60	Friedel–Crafts Acylation	537
61	The Analysis of Antihistamine Drugs by Gas Chromatography–Mass Spectrometry	545

- 62 The Use of Organozinc Reagents in Synthesis: An Exercise in Synthesis and Structure Proof by Spectroscopy 548
- 63 Synthesis of Naproxen by Palladium Catalysis 552
- 64 The Aldehyde Engima 566
- 65 Synthesis of Substituted Chalcones: A Guided-Inquiry Experience 569
- 66 Green Epoxidation of Chalcones 574
- 67 Cyclopropanation of Chalcones 578
- 68 Michael and Aldol Condensation Reactions 582
- 69 Esterification Reactions of Vanillin: The Use of NMR to Solve a Structure Proof Problem 586

PART 6 The Techniques 589

- 1 Laboratory Safety 590
- 2 The Laboratory Notebook, Calculations, and Laboratory Records 609
- 3 Laboratory Glassware: Care and Cleaning 617
- 4 How to Find Data for Compounds: Handbooks and Catalogs 625
- 5 Measurement of Volume and Weight 632
- 6 Heating and Cooling Methods 640
- 7 Reaction Methods 647
- 8 Filtration 667
- 9 Physical Constants of Solids: The Melting Point 678
- 10 Solubility 687
- 11 Crystallization: Purification of Solids 696
- 12 Extractions, Separations, and Drying Agents 718
- 13 Physical Constants of Liquids: The Boiling Point and Density 745
- 14 Simple Distillation 756
- 15 Fractional Distillation, Azeotropes 768
- 16 Vacuum Distillation, Manometers 785
- 17 Sublimation 797
- 18 Steam Distillation 802
- 19 Column Chromatography 808
- 20 Thin-Layer Chromatography 828
- 21 High-Performance Liquid Chromatography (HPLC) 842
- 22 Gas Chromatography 847
- 23 Polarimetry 867
- 24 Refractometry 875
- 25 Infrared Spectroscopy 880
- 26 Nuclear Magnetic Resonance Spectroscopy (Proton NMR) 914
- 27 Carbon-13 Nuclear Magnetic Resonance Spectroscopy 921
- 28 Mass Spectrometry 969
- 29 Guide to the Chemical Literature 987

Appendices 1001

- 1 Tables of Unknowns and Derivatives 1002
- 2 Procedures for Preparing Derivatives 1016
- 3 Index of Spectra 1020

Index 1023

Introduction to Basic Laboratory Techniques

Introduction to Microscale Laboratory

This textbook discusses the important laboratory techniques of organic chemistry and illustrates many important reactions and concepts. In the traditional approach to teaching this subject, the quantities of chemicals used were on the order of 5–100 grams, and glassware was designed to contain up to 500 mL of liquid. This scale of experiment we might call a **macroscale** experiment. The approach used here, a **microscale** approach, differs from the traditional laboratory course in that nearly all the experiments use small amounts of chemicals. Quantities of chemicals used range from about 50 to 1000 *milligrams* (0.050–1.000 g), and glassware is designed to contain less than 25 mL of liquid. The advantages include improved safety in the laboratory, reduced risk of fire and explosion, and reduced exposure to hazardous vapors. This approach decreases the need for hazardous waste disposal, leading to reduced contamination of the environment. You will learn to work with the same level of care and neatness that has previously been confined to courses in analytical chemistry.

This experiment introduces the equipment and shows how to construct some of the apparatus needed to carry out further experiments. Detailed discussion of how to assemble apparatus and how to practice the techniques is found in Part Six (“The Techniques”) of this textbook. This experiment provides only a brief introduction, sufficient to allow you to begin working. You will need to read the techniques chapters for more complete discussions.

Microscale organic experiments require you to develop careful laboratory techniques and to become familiar with apparatus that is somewhat unusual, compared with traditional glassware. We strongly recommend that each student do Laboratory Exercises 1 and 2. These exercises will acquaint you with the most basic microscale techniques. To provide a strong foundation, we further recommend that each student complete most of Experiments 2 through 18 in Part One of this textbook before attempting any other experiments in the textbook.

READ Technique 1 “Laboratory Safety.”

HEATING METHODS

Aluminum Block

The most convenient means of heating chemical reactions on a small scale is to use an **aluminum block**. An aluminum block consists of a square of aluminum that has holes drilled into it. The holes are sized to correspond to the diameters of the most common vials and flasks that are likely to be heated. Often there is also

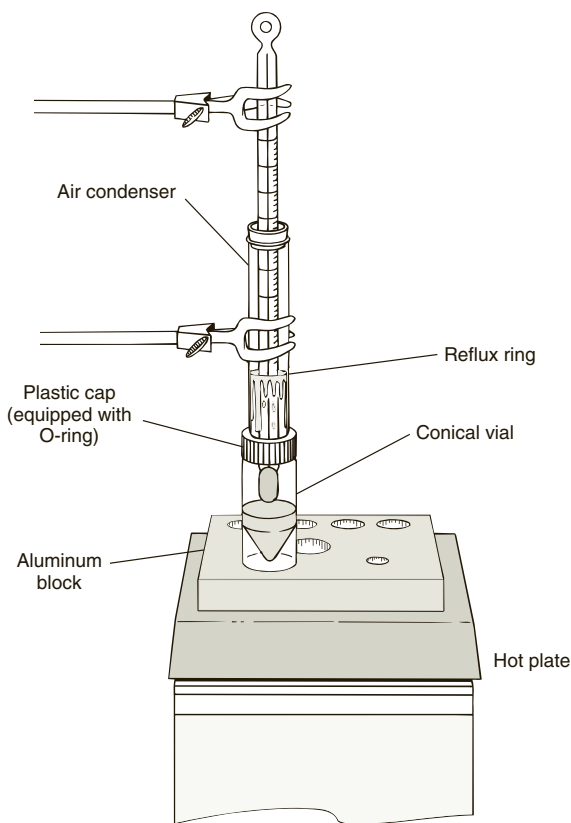


Figure 1
Aluminum block with hot plate and reflux apparatus.

a hole intended to accept the bulb of a thermometer, so that the temperature of the block can be monitored. However, this practice is not recommended. The aluminum block is heated by placing it on a hot plate. An aluminum block is shown in Figure 1. Note that the thermometer in this figure is not used to monitor the temperature of the block.

CAUTION



You should not use a mercury thermometer in direct contact with an aluminum block. If it breaks, the mercury will vaporize on the hot surface. Instead, use a nonmercury thermometer, a metal dial thermometer, or a digital electronic temperature-measuring device. See Technique 6, Section 6.1.

It is recommended that an equipment kit contain two aluminum blocks, one drilled with small holes and able to accept the conical vials found in the glassware kit and another drilled with larger holes and able to accept small round-bottom flasks. The aluminum blocks can be made from inexpensive materials in a small mechanical shop, or they can be purchased from a glassware supplier.

Sand Baths

Another commonly used means of heating chemical reactions on a small scale is to use a **sand bath**. The sand bath consists of a Petri dish or a small crystallizing dish that has been filled to a depth of about 1 cm with sand. The sand bath is also heated by placing it on a hot plate. The temperature of the sand bath may be monitored by clamping a thermometer in position so that the bulb of the thermometer is buried in the sand. A sand bath, with thermometer, is shown in Figure 2.

We recommend that an aluminum block, rather than a sand bath, be used as a heating source whenever possible. The aluminum block can be heated and cooled quickly, it is indestructible, and there are no problems with spillage of sand.

Water Bath

When precise control at lower temperatures (below about 80°C) is desired, a suitable alternative is to prepare a **water bath**. The water bath consists of a beaker filled to the required depth with water. The hot plate is used to heat the water bath to the desired temperature. The water in the water bath can evaporate during heating. It is useful to cover the top of the beaker with aluminum foil to diminish this problem.

CONICAL REACTION VIALS

One of the most versatile pieces of glassware contained in the microscale organic glassware kit is the **conical reaction vial**. This vial is used as a vessel in which organic reactions are performed. It may serve as a storage container. It is also used for extractions (see Technique 12). A reaction vial is shown in Figure 3.

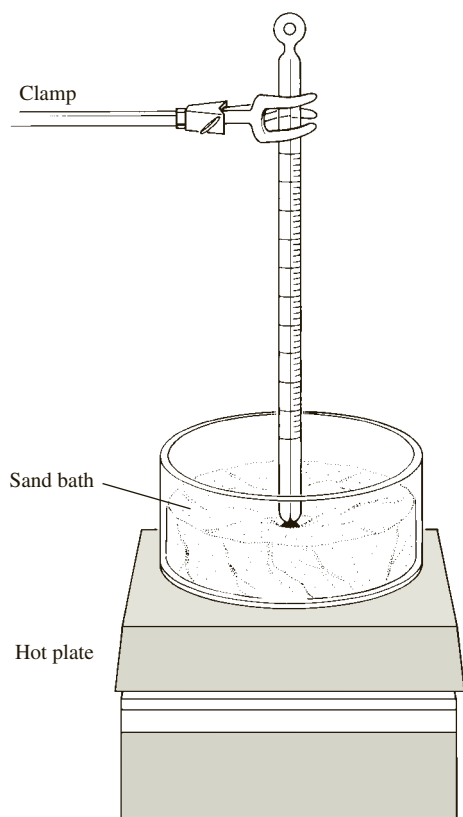


Figure 2
Sand bath with hot plate and thermometer.

The flat base of the vial allows it to stand upright on the laboratory bench. The interior of the vial tapers to a narrow bottom. This shape makes it possible to withdraw liquids completely from the vial, using a disposable Pasteur pipette. The vial has a screw cap, which tightens by means of threads cast into the top of the vial. The top also has a ground-glass inner surface. This ground-glass joint allows you to assemble components of glassware tightly.

The plastic cap that fits the top of the conical vial has a hole in the top. This hole is large enough to permit the cap to fit over the inner joints of other components of the glassware kit (see Figure 4). A Teflon insert, or liner, fits inside the cap to cover the hole when the cap is used to seal a vial tightly. Notice that only one side of the liner is coated with Teflon; the other side is coated with a silicone rubber. The Teflon side generally is the harder side of the insert, and it will feel more slippery. The Teflon side should always face toward the inside of the vial. An O-ring fits inside the cap when the cap is used to fasten pieces of glassware together. The cap and its Teflon insert are shown in the expanded view in Figure 3.

NOTE: Do not use the O-ring when the cap is used to seal the vial.

You can assemble the components of the glassware kit into one unit that holds together firmly and clamps easily to a ring stand. Slip the cap from the conical vial over the inner (male) joint of the upper piece of glassware and fit a rubber O-ring over the inner joint. Then assemble the apparatus by fitting the inner ground-glass joint into the outer (female) joint of the reaction vial and tighten the screw cap to attach the entire apparatus firmly together. The assembly is illustrated in Figure 4.

The walls of the conical vials are made of thick glass. Heat does not transfer through these walls very quickly. This means that if the vial is subjected to rapid changes in temperature, strain building up within the glass walls of the vial may cause the glass to crack. For this reason, do not attempt to cool these vials quickly by running cold water on them. It is safer to allow them to cool naturally by allowing them to stand.

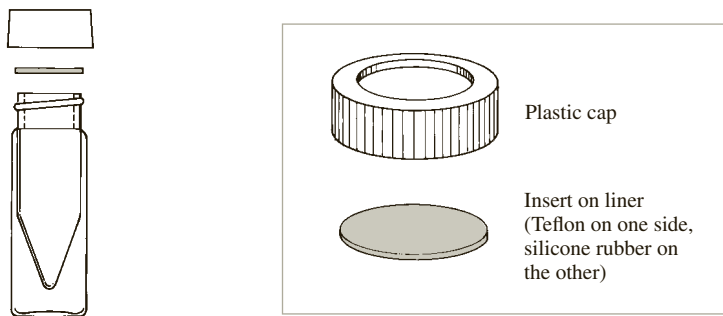


Figure 3
A conical reaction vial. (The inset shows an expanded view of the cap with its Teflon insert.)

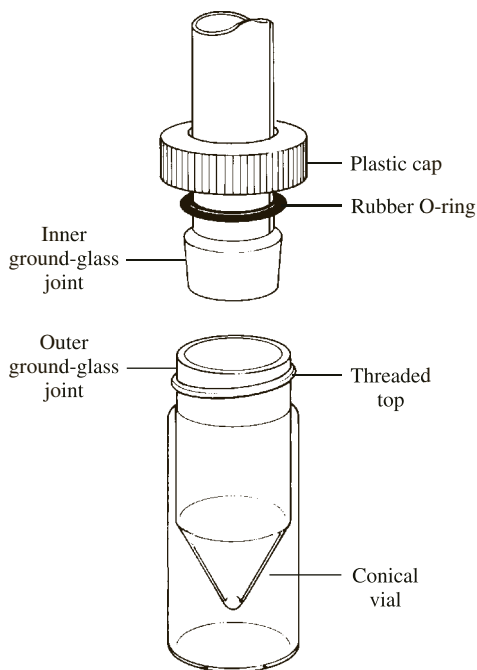


Figure 4
Assembling glassware components.

Although the conical vials have flat bottoms intended to allow them to stand up on the laboratory bench, this does not prevent them from falling over.

NOTE: It is good practice to store the vials standing upright inside small beakers.

The vials are somewhat top-heavy, and it is easy to upset them. The beaker will prevent the vial from falling over onto its side.

MEASUREMENT OF SOLIDS

Weighing substances to the nearest milligram requires that the weighings be done on a sensitive **top-loading balance** or an **analytical balance**.

NOTE: You must not weigh chemicals directly on balance pans.

Many chemicals can react with the metal surface of the balance pan and thus ruin it. All weighings must be made into a container that has been weighed previously (**tared**). This tare weight is subtracted from the total weight of container plus sample to give the weight of the sample. Some balances have a built-in compensating feature, the tare button, that allows you to subtract the tare weight of the container automatically, thus giving the weight of the sample directly. A top-loading and an analytical balance are shown in Figure 5.

Balances of this type are quite sensitive and expensive. Take care not to spill chemicals on the balance. It is also important to make certain that any spilled materials are cleaned up immediately.

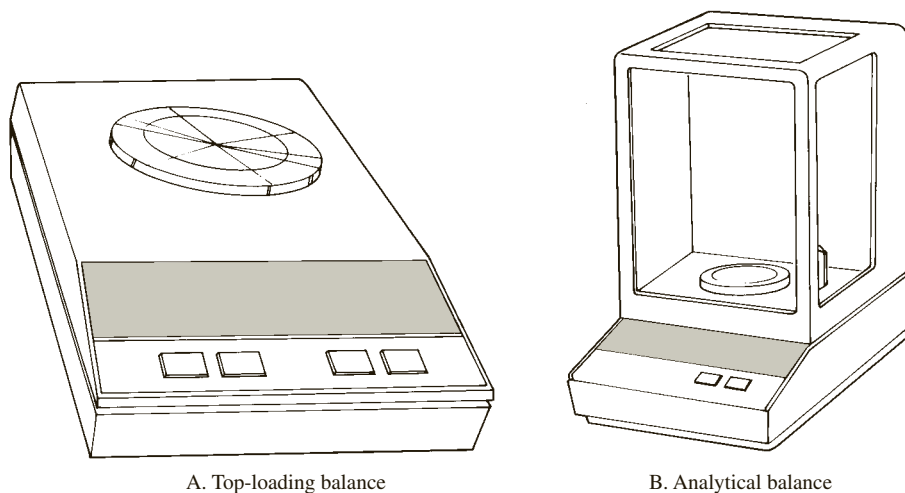


Figure 5
Laboratory balances.

MEASUREMENT OF LIQUIDS

In microscale experiments, liquid samples are measured using a pipette. When small quantities are used, graduated cylinders do not provide the accuracy needed to give good results. There are two common methods of delivering known amounts of liquid samples, **automatic pipettes** and **graduated pipettes**. When accurate quantities of liquid reagents are required, the best technique is to deliver the desired amount of liquid reagent from the pipette into a container whose tare weight has been determined previously. The container, with sample, is then weighed a second time in order to obtain a precise value of the amount of reagent.

Automatic Pipettes

Automatic pipettes may vary in design, according to the manufacturer. The following description, however, should apply to most models. The automatic pipette consists of a handle that contains a spring-loaded plunger and a micrometer dial. The dial controls the travel of the plunger and is the means used to select the amount of liquid that the pipette is intended to dispense. Automatic pipettes are designed to deliver liquids within a particular range of volumes. For example, a pipette may be designed to cover the range from 10 to 100 μL (0.010 to 0.100 mL) or from 100 to 1000 μL (0.100 to 1.000 mL).

Automatic pipettes must never be dipped directly into the liquid sample without a plastic tip. The pipette is designed so that the liquid is drawn only into the tip. The liquids are never allowed to come in contact with the internal parts of the pipette. The plunger has two **detent**, or "stop," positions used to control the filling and dispensing steps. Most automatic pipettes have a stiffer spring that controls the movement of the plunger from the first to the second detent position. You will find a greater resistance as you press the plunger past the first detent.

To use the automatic pipette, follow the steps as outlined here. These steps are also illustrated in Figure 6.

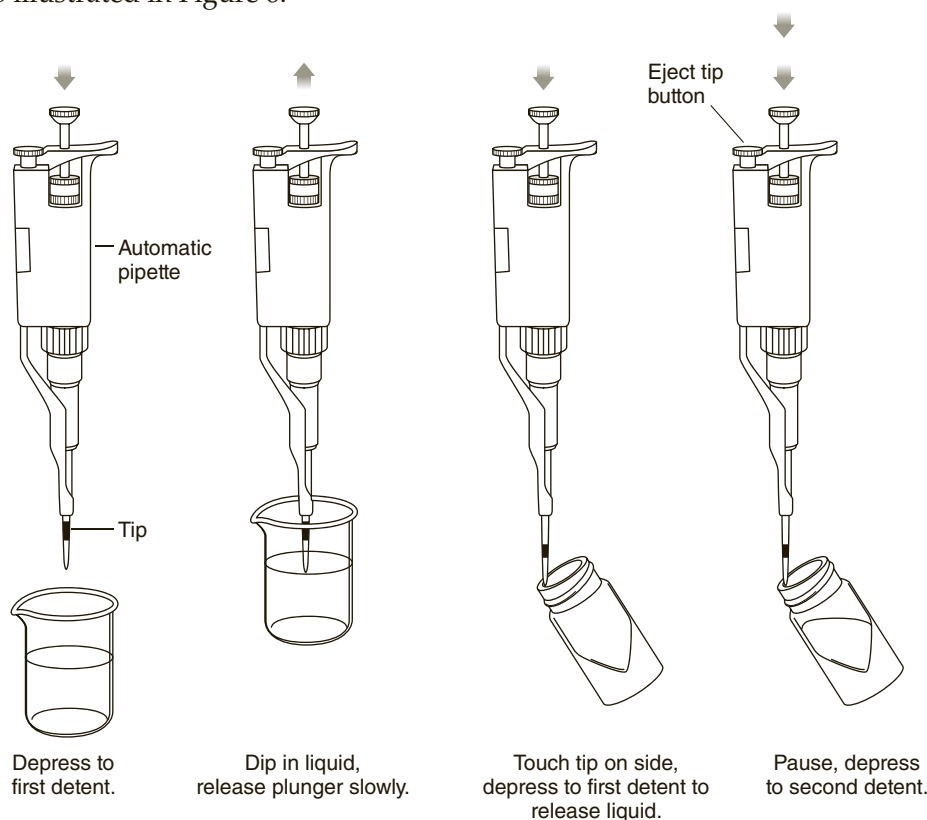


Figure 6
Use of an automatic pipette.

1. Select the desired volume by adjusting the micrometer control on the pipette handle.
2. Place a plastic tip on the pipette. Be certain that the tip is attached securely.
3. Push the plunger down to the first detent position. Do not press the plunger to the second position. If you press the plunger to the second detent, an incorrect volume of liquid will be delivered.
4. Dip the tip of the pipette into the liquid sample. Do not immerse the entire length of the plastic tip in the liquid. It is best to dip the tip only to a depth of about 1 cm.
5. Release the plunger *slowly*. Do not allow the plunger to snap back, or liquid may splash up into the plunger mechanism and ruin the pipette. Furthermore, rapid release of the plunger may cause air bubbles to be drawn into the pipette. At this point, the pipette has been filled.
6. Move the pipette to the receiving vessel. Touch the tip of the pipette to an interior wall of the container.
7. Slowly push the plunger down to the first detent. This action dispenses the liquid into the container.
8. Pause 1–2 seconds and then depress the plunger to its second detent position to expel the last drop of liquid. The action of the plunger may be stiffer in this range than it was up to the first detent.
9. Withdraw the pipette from the receiver. If the pipette is to be used with a different liquid, remove the pipette tip and discard it.

Automatic pipettes are designed to deliver aqueous solutions with an accuracy of within a few percentage points. The amount of liquid actually dispensed varies, however, depending on the viscosity, surface tension, and vapor pressure of the liquid. The typical automatic pipette is very accurate with aqueous solutions but is not always as accurate with other liquids.

Dispensing Pumps

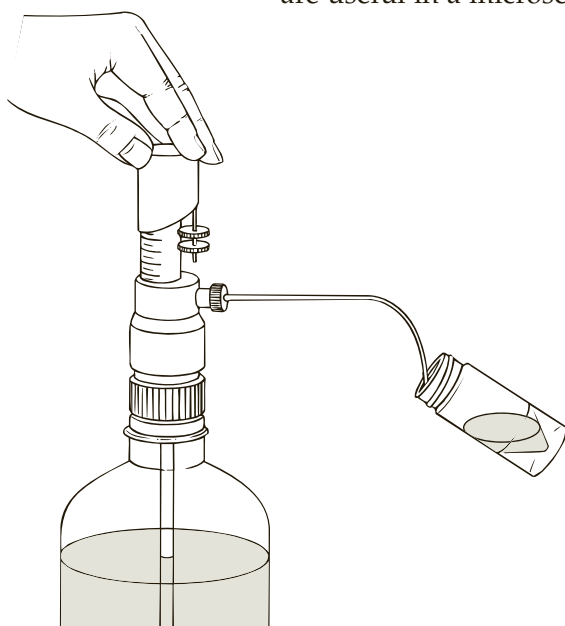


Figure 7
Use of a dispensing pump.

Some scientific supply catalogs offer a series of dispensing pumps. These pumps are useful in a microscale organic laboratory because they are simple to operate, easy to clean, chemically inert, and quite accurate. The interior parts of dispensing pumps are made of Teflon, which renders them inert to most organic solvents and reagents. A dispensing pump is illustrated in Figure 7.

The first step in using a dispensing pump is to adjust the pump so that it dispenses the desired volume of liquid. Normally, the instructor will make this adjustment. Once the pump is adjusted correctly, it is a simple matter to dispense a liquid. Simply lift the head of the pump as far as it will travel. When you release the head, it will fall, and the liquid will issue from the spout. With viscous liquids, the head of the pump may not fall by itself. In such an instance, gently guide the head downward. After the liquid has been dispensed, you should touch the tip of the dispensing tube to an interior wall of the container in order to remove the last drop of liquid.

As with automatic pipettes, dispensing pumps are designed to deliver aqueous solutions with an accuracy of within a few percentage points. The amount of liquid actually dispensed will vary, however, depending on

the viscosity, surface tension, and vapor pressure of the liquid. You should always weigh the liquid to determine the amount accurately.

Graduated Pipettes

A less-expensive means of delivering known quantities of liquid is to use a graduated pipette. Graduated pipettes should be familiar to those of you who have taken general chemistry or quantitative analysis courses. Because they are made of glass, they are inert to most organic solvents and reagents. Disposable serological pipettes may be an attractive alternative to standard graduated pipettes. The 2-mL size of a disposable pipette represents a convenient size for the organic laboratory.

Never draw liquids into the pipettes using mouth suction. A pipette bulb or a pipette pump, not a rubber dropper bulb, must be used to fill pipettes. We recommend the use of a pipette pump. A pipette fits snugly into the pipette pump, and the pump can be controlled to deliver precise volumes of liquids. Control of the pipette pump is accomplished by rotating a knob on the pump. Suction created when the knob is turned draws the liquid into the pipette. Liquid is expelled from the pipette by turning the knob in the opposite direction. The pump works satisfactorily with organic, as well as aqueous, solutions.

An alternative, and less expensive, approach is to use a rubber pipette bulb. Use of the pipette bulb is made more convenient by inserting a plastic automatic pipette tip into a rubber pipette bulb.¹ The tapered end of the pipette tip fits snugly into the end of a pipette. Drawing the liquid into the pipette is made easy, and it is also convenient to remove the pipette bulb and place a finger over the pipette opening to control the flow of liquid.

The calibrations printed on graduated pipettes are reasonably accurate, but you should practice using the pipettes in order to achieve this accuracy. When accurate quantities of liquids are required, the best technique is to weigh the reagent that has been delivered from the pipette.

LABORATORY EXERCISE 1

Option A, Automatic Pipette

Accurately weigh a 3-mL conical vial, with screw cap and Teflon insert, on a balance. Determine its weight to the nearest milligram (nearest 0.001 g). Using the automatic pipette, dispense 0.500 mL of water into the vial, replace the cap assembly (with the insert arranged Teflon side down), and weigh the vial a second time. Determine the weight of water dispensed. Calculate the density of water from your results. Repeat the experiment using 0.500 mL of hexane. Dispose of any excess hexane in a designated waste container. Calculate the density of hexane from your data. Record the results in your notebook, along with your comments on any deviations from literature values that you may have noticed. At room temperature, the density of water is 0.997 g/mL, and the density of hexane is 0.660 g/mL.

Option B, Dispensing Pump

Accurately weigh a 3-mL conical vial, with screw cap and Teflon insert, on a balance. Determine its weight to the nearest milligram (nearest 0.001 g). Using a dispensing pump that has been adjusted to deliver 0.500 mL, dispense 0.500 mL of water into the vial, replace the cap assembly, and weigh the vial a second time. Determine the weight of water dispensed. Calculate the density of water from your results. Repeat the experiment using 0.500 mL of hexane. Dispose of any excess hexane in a designated waste container. Calculate the density of hexane from your data.

¹ This technique was described in Deckey, G. A Versatile and Inexpensive Pipet Bulb. *Journal of Chemical Education*, 57 (July 1980): 526.

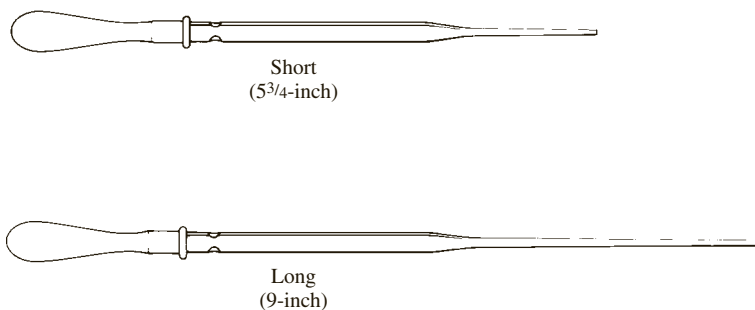


Figure 8
Disposable Pasteur pipettes.

Record the results in your notebook, along with your comments on any deviations from literature values that you may have noticed. See Option A for the density of water and of hexane.

Option C, Graduated Pipette

Accurately weigh a 3-mL conical vial, with screw cap and Teflon insert, on a balance. Determine its weight to the nearest milligram (nearest 0.001 g). Using a 1.0-mL graduated pipette, dispense 0.50 mL of water into the vial, replace the cap assembly, and weigh the vial a second time. Determine the weight of water dispensed. Calculate the density of water from your results. Repeat the experiment using 0.50 mL of hexane. Dispose of any excess hexane in a designated waste container. Calculate the density of hexane from your data. Record the results in your notebook along with your comments on any deviations from literature values that you may have noticed. See Option A for the density of water and of hexane

Disposable (Pasteur) Pipettes

A convenient way of dispensing liquids when a great deal of accuracy is not required is to use a disposable pipette, or Pasteur pipette. Two sizes of Pasteur pipettes are shown in Figure 8. Even though accurate calibration may not be required when these pipettes are used, it is nevertheless handy to have some idea of the volume contained in the pipette. A crude calibration is, therefore, recommended.

LABORATORY EXERCISE 2

Pipette Calibration

On a balance, weigh 0.5 grams (0.5 mL) of water into a 3-mL conical vial. Select a short (5 $\frac{3}{4}$ -inch) Pasteur pipette and attach a rubber bulb. Squeeze the rubber bulb before inserting the tip of the pipette into the water. Try to control how much you depress the bulb, so that when the pipette is placed into the water and the bulb is completely released, only the desired amount of liquid is drawn into the pipette. (This skill may take some time to acquire, but it will facilitate your use of a Pasteur pipette.) When the water has been drawn up, place a mark with an indelible marking pen at the position of the meniscus. A more durable mark can be made by scoring the pipette with a file. Repeat this procedure with 1.0 gram of water, and make a 1-mL mark on the same pipette.

Additional Pasteur pipettes can be calibrated easily by holding them next to the pipette calibrated in Laboratory Exercise 2 and scoring a new mark on each pipette at the same level as the mark placed on the calibrated pipette. We recommend that several Pasteur pipettes be calibrated at one time for use in future experiments.

Extraction

A technique frequently applied in purifying organic reaction products is **extraction**. In this method, a solution is mixed thoroughly with a second solvent. The second

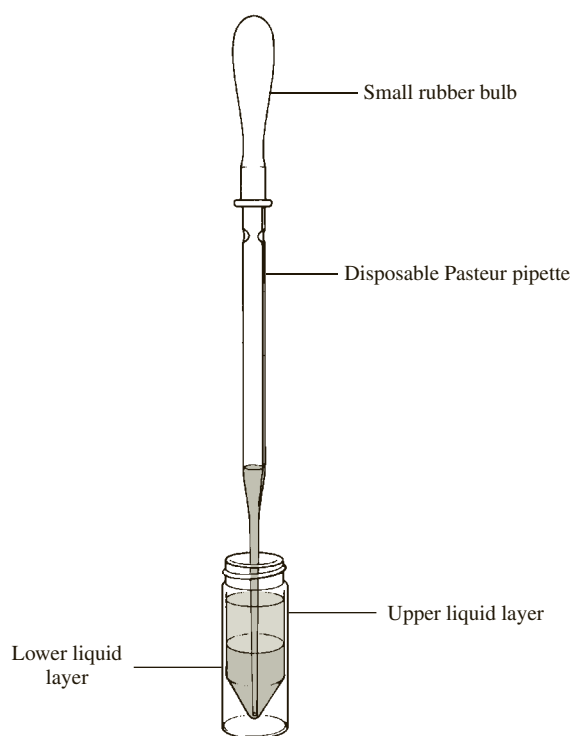


Figure 9
Separation of immiscible liquid layers in a conical vial.

solvent is not miscible with the first solvent. When the two solvents are mixed, the dissolved substances (solutes) distribute themselves between the two solvents until an equilibrium is established. When the mixing is stopped, the two immiscible solvents separate into two distinct layers. The solutes are distributed between the two solvents so that each solute is found in greater concentration in that solvent in which it is more soluble. Separation of the two immiscible solvent layers thus becomes a means of separating solutes from one another based on their relative solubilities in the two solvents.

In a common application, an aqueous solution may contain both inorganic and organic products. An organic solvent that is immiscible with water is added, and the mixture is shaken thoroughly. When the two solvent layers are allowed to form again, on standing, the organic solutes are transferred to the organic solvent while the inorganic solutes remain in the aqueous layer. When the two layers separate, the organic and inorganic products also separate from one another. The separation, as described here, may not be complete. The inorganic materials may be somewhat soluble in the organic solvent, and the organic products may retain some water solubility. Nevertheless, reasonably complete separations of reaction products can be achieved by the extraction method.

For microscale experiments, the conical reaction vial is the glassware item used for extractions. Place the two immiscible liquid layers in the vial, and seal the top with a screw cap and a Teflon insert (Teflon side toward the inside of the vial). Shake the vial to provide thorough mixing between the two liquid phases. As the shaking continues, vent the vial periodically by loosening the cap and then tightening it again. After about 5 or 10 seconds of shaking, loosen the cap to vent the vial, retighten it, and allow the vial to stand upright in a beaker until the two liquid layers separate completely.

The two liquid layers are separated by withdrawing the *lower* layer using a disposable Pasteur pipette. This separation technique is illustrated in Figure 9. Take care not to disturb the liquid layers by allowing bubbles to issue from the pipette. Squeeze the pipette bulb to the required amount before introducing the pipette into the vial. Also take care not to allow any of the upper liquid layer to enter the pipette. The pointed shape of the interior of the conical vial makes it easy to remove all the lower layer without allowing it to be contaminated by some of the upper liquid layer. More precise control in the separation can be achieved by using a filter-tip pipette (see Technique 8, Section 8.6).

Other Useful Techniques

The practice of organic chemistry requires you to master many more techniques than the ones described in this experiment. Those techniques included here are only the most elementary ones, those needed to get you started in the laboratory. Additional techniques are described fully in Part Six of this textbook, and Experiments 2 through 18 expose you to the most important of them.

Some other practical hints need to be introduced at this point. The first of these involves manipulating small amounts of solid substances. The efficient transfer of solids requires a small spatula. We recommend that you have two **microspatulas**, similar to those shown in Figure 10, as part of your standard desk stock. The design

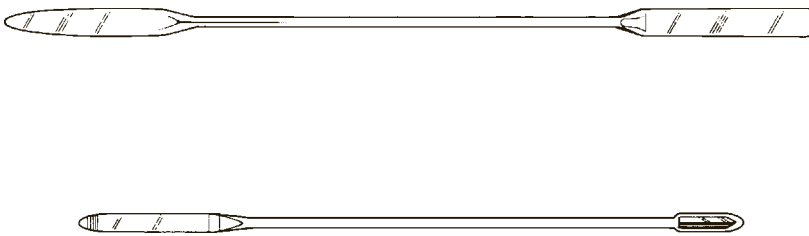


Figure 10
Microspatulas.

of these spatulas permits the handling of milligram quantities of substances without undue spillage or waste. The larger style (see Figure 10) is more useful when relatively large quantities of solid must be dispensed.

A clean work area is of utmost importance when working in the laboratory. The need for cleanliness is particularly great when working with the small amounts of materials characteristic of microscale laboratory experiments.

NOTE: You must read Technique 1 “Laboratory Safety.” In preventing accidents, there is no substitute for preparation and care.

With this final word of caution and advice, we hope you enjoy the learning experience you are about to begin. Learning the care and precision that microscale experiments require may seem difficult at first, but before long you will be comfortable with the scale of the experiments. You will develop much better laboratory technique as a result of microscale practice, and this added skill will serve you well.

Solubility

Solubility

Polarity

Acid-base chemistry

Critical thinking application

Having a good comprehension of solubility behavior is essential for understanding many procedures and techniques in the organic chemistry laboratory. For a thorough discussion of solubility, read the chapter on this concept (Technique 10) before proceeding because an understanding of this material is assumed in this experiment.

In Parts A and B of this experiment, you will investigate the solubility of various substances in different solvents. As you are performing these tests, it is helpful to pay attention to the polarities of the solutes and solvents and to even make predictions based on this (see “Guidelines for Predicting Polarity and Solubility,” Technique 10, Section 10.4). The goal of Part C is similar to that of Parts A and B, except that you will be looking at miscible and immiscible pairs of liquids. In Part D, you will investigate the solubility of organic acids and bases. Technique 10, Section 10.2B will help you understand and explain these results.

In Part E, you will perform several exercises that involve the application of the solubility principles learned in Parts A–D of this experiment.

REQUIRED READING

New: Technique 5 Measurement of Volume and Weight
Technique 10 Solubility

SUGGESTED WASTE DISPOSAL

Dispose of all wastes containing methylene chloride into the container marked for halogenated waste. Place all other organic wastes into the nonhalogenated organic waste container.

NOTES TO THE INSTRUCTOR

In Part A of the procedure, it is important that students follow the instructions carefully. Otherwise, the results may be difficult to interpret. It is particularly important

that consistent stirring be done for each solubility test. This can be done most easily by using the larger-style microspatula shown in Figure 10 in Experiment 1, Introduction to Microscale Laboratory.

We have found that some students have difficulty performing Critical Thinking Application 2 in Part E on the same day that they complete the rest of this experiment. Many students need time to assimilate the material in this experiment before they can complete this exercise successfully. One approach is to assign Critical Thinking Applications from several technique experiments (for example, Experiments 2–4) to a laboratory period after students complete the individual technique experiments. This provides an effective way of reviewing some of the basic techniques.

PROCEDURE

NOTE: It is very important that you follow these instructions carefully and that consistent stirring be done for each solubility test.




Part A. Solubility of Solid Compounds

Place about 40 mg (0.040 g) of benzophenone into each of four *dry* test tubes.¹ (Don't try to be exact: You can be 1–2 mg off and the experiment will still work.) Label the test tubes and then add 1 mL of water to the first tube, 1 mL of methyl alcohol to the second tube, and 1 mL of hexane to the third tube. The fourth tube will serve as a control. Determine the solubility of each sample in the following way: Using the rounded end of a microspatula (the larger style, in Figure 10 in Experiment 1, Introduction to Microscale Laboratory), stir each sample continuously for 60 seconds by twirling the spatula rapidly. If a solid dissolves completely, note how long it takes for the solid to dissolve. *After 60 seconds* (do not stir longer), note whether the compound is soluble (dissolves completely), insoluble (none of it dissolves), or partially soluble. You should compare each tube with the control in making these determinations. You should state that a sample is partially soluble only if a significant amount (at least 50%) of the solid has dissolved. If it is not clear that a significant amount of solid has dissolved, then state that the sample is insoluble. If all but a couple of granules have dissolved, state that the sample is soluble. An additional hint for determining partial solubility is given in the next paragraph. Record these results in your notebook in the form of a table, as shown on this page. For those substances that dissolve completely, note how long it took for the solid to dissolve.

Although the instructions just given should enable you to determine whether a substance is partially soluble, you may use the following procedure to confirm this. Using a Pasteur pipette, carefully remove most of the solvent from the test tube *while leaving the solid behind*. Transfer the liquid to another test tube and then evaporate the solvent by heating the tube in a hot water bath. Directing a stream of air or nitrogen gas into the tube will speed up the evaporation (see Technique 7, Section 7.10). When the solvent has completely evaporated, examine the test tube for any remaining solid. If there is solid in the test tube, the compound is partially soluble. If there is no, or very little, solid remaining, you can assume that the compound is insoluble.

Now repeat the directions just given, substituting malonic acid and biphenyl for benzophenone. Record these results in your notebook.

¹ *Note to the instructor:* Grind up the benzophenone flakes into a powder.

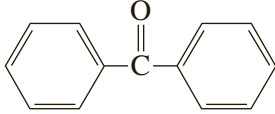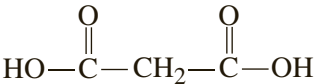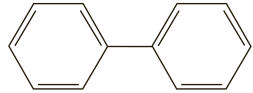Organic Compounds	Solvents		
	Water (highly polar)	Methyl Alcohol (intermediate polarity)	Hexane (nonpolar)
Benzophenone 			
Malic acid 			
Biphenyl 			

Part B. Solubility of Different Alcohols

For each solubility test (see table below), add 1 mL of solvent (water or hexane) to a test tube. Then add one of the alcohols, *dropwise*. Carefully observe what happens as you add each drop. If the liquid solute is soluble in the solvent, you may see tiny horizontal lines in the solvent. These mixing lines indicate that solution is taking place. *Shake the tube after adding each drop*. While you shake the tube, the liquid that was added may break up into small balls that disappear in a few seconds. This also indicates that solution is taking place. Continue adding the alcohol with shaking until you have added a total of 20 drops. If an alcohol is partially soluble, you will observe that at first the drops will dissolve, but eventually a second layer of liquid (undissolved alcohol) will form in the test tube. Record your results (soluble, insoluble, or partially soluble) in your notebook in table form.

Alcohols	Solvents	
	Water	Hexane
1-Octanol $\text{CH}_3(\text{CH}_2)_6\text{CH}_2\text{OH}$		
1-Butanol $\text{CH}_3\text{CH}_2\text{CH}_2\text{CH}_2\text{OH}$		
Methyl alcohol CH_3OH		

Part C. Miscible or Immiscible Pairs

For each of the following pairs of compounds, add 1 mL of each liquid to the same test tube. Use a different test tube for each pair. Shake the test tube for 10–20 seconds to determine whether the two liquids are miscible (form one layer) or immiscible (form two layers). Record your results in your notebook.

Water and ethyl alcohol

Water and diethyl ether

Water and methylene chloride



Water and hexane

Hexane and methylene chloride

Part D. Solubility of Organic Acids and Bases

Place about 30 mg (0.030 g) of benzoic acid into each of three *dry* test tubes. Label the test tubes and then add 1 mL of water to the first tube, 1 mL of 1.0 M NaOH to the second tube, and 1 mL of 1.0 M HCl to the third tube. Stir the mixture in each test tube with a microspatula for 10–20 seconds. Note whether the compound is soluble (dissolves completely) or is insoluble (none of it dissolves). Record these results in the table. Now take the second tube containing benzoic acid and 1.0 M NaOH. With stirring add 6 M HCl dropwise until the mixture is acidic. Test the mixture with litmus or pH paper to determine when it is acidic.² When it is acidic, stir the mixture for 10–20 seconds and note the result (soluble or insoluble) in the table.

Repeat this experiment using ethyl 4-aminobenzoate and the same three solvents. Record the results. Now take the third tube containing ethyl 4-aminobenzoate and 1.0 M HCl. With stirring, add 6 M NaOH dropwise until the mixture is basic. Test the mixture with litmus or pH paper to determine when it is basic. Stir the mixture for 10–20 seconds and note the result.

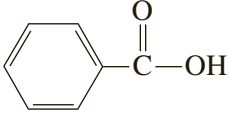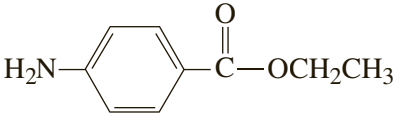Compounds	Solvents		
	Water	1.0 M NaOH	1.0 M HCl
Benzoic acid 		Add 6.0 M HCl	
Ethyl 4-aminobenzoate 			Add 6.0 M NaOH

Part E. Critical Thinking Applications

- Determine by experiment whether each of the following pairs of liquids is miscible or immiscible.

Acetone and water

Acetone and hexane


How can you explain these results, given that water and hexane are immiscible?

- You will be given a test tube containing two immiscible liquids and a solid organic compound that is dissolved in one of the liquids.³ You will be told the identity of the two liquids and the solid compound, but you will not know the relative positions of the two liquids or in which liquid the solid is dissolved.

² Do not place the litmus or pH paper into the sample; the dye will dissolve. Instead, place a drop of solution from your spatula onto the test paper. With this method, several tests can be performed using a single strip of paper.

³ The sample you are given may contain one of the following combinations of solid and liquids (the solid is listed first): fluorene, methylene chloride, water; triphenylmethanol, diethyl ether, water; salicylic acid, methylene chloride, 1 M NaOH; ethyl 4-aminobenzoate, diethyl ether, 1 M HCl; naphthalene, hexane, water; benzoic acid, diethyl ether, 1 M NaOH; *p*-aminoacetophenone, methylene chloride, 1 M HCl.

Consider the following example, in which the liquids are water and hexane and the solid compound is biphenyl.

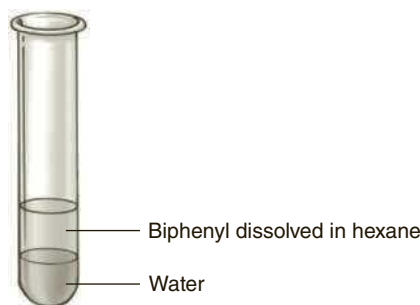


- a. Without doing any experimental work, predict where each liquid is (top or bottom) and in which liquid the solid is dissolved. Justify your prediction. You may want to consult a handbook such as *The Merck Index* or the *CRC Handbook of Chemistry and Physics* to determine the molecular structure of a compound or to find any other relevant information. Note that dilute solutions such as 1 M HCl are composed mainly of water, and the density will be close to 1.0 g/mL. Furthermore, you should assume that the density of a solvent is not altered significantly when a solid dissolves in the solvent.
 - b. Now try to prove your prediction experimentally. That is, demonstrate which liquid the solid compound is dissolved in and the relative positions of the two liquids. You may use any experimental technique discussed in this experiment or any other technique that your instructor will let you try. In order to perform this part of the experiment, it may be helpful to separate the two layers in the test tube. This can be done easily and effectively with a Pasteur pipette. Squeeze the bulb on the Pasteur pipette and then place the tip of the pipette on the bottom of the test tube. Now withdraw only the bottom layer and transfer it to another test tube. Note that evaporating the water from an aqueous sample takes a very long time; therefore, this may not be a good way to show that an aqueous solution contains a dissolved compound. However, other solvents may be evaporated more easily. Explain what you did and whether or not the results of your experimental work were consistent with your prediction.
3. Add 0.025 g of tetraphenylcyclopentadienone to a dry test tube. Add 1 mL of methyl alcohol to the tube and shake for 60 seconds. Is the solid soluble, partially soluble, or insoluble? Explain your answer based on the experimental results. Repeat this test using 1-butanol, rather than methyl alcohol as the solvent. Explain the different results.

REPORT

Part A

1. Summarize your results in table form.
2. Explain the results for all the tests done. In explaining these results, you should consider the polarities of the compound and the solvent and the potential for hydrogen bonding. For example, consider a similar solubility test for *p*-dichlorobenzene in hexane. The test indicates that *p*-dichlorobenzene is soluble in hexane. This result can be explained by stating that hexane is nonpolar,

whereas *p*-dichlorobenzene is slightly polar. Because the polarities of the solvent and solute are similar, the solid is soluble. (Remember that the presence of a halogen does not significantly increase the polarity of a compound.)

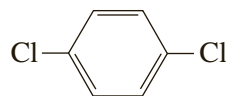

p-Dichlorobenzene

- There should be a difference in your results between the solubilities of biphenyl and benzophenone in methyl alcohol. Explain this difference.
- There should be a difference in your results between the solubilities of benzophenone in methyl alcohol and benzophenone in hexane. Explain this difference.

Part B

- Summarize your results in table form.
- Explain the results for the tests done in water. In explaining these results, you should consider the polarities of the alcohols and water.
- Explain, in terms of polarities, the results for the tests done in hexane.

Part C

- Summarize your results in table form.
- Explain the results in terms of polarities and/or hydrogen bonding.

Part D


- Summarize your results in table form.
- Explain the results for the tube in which 1.0 M NaOH was added to benzoic acid. Write an equation for this, using complete structures for all organic substances. Now describe what happened when 6.0 M HCl was added to this same tube and explain this result.
- Explain the results for the tube in which 1.0 M HCl was added to ethyl 4-aminobenzoate. Write an equation for this. Now describe what happened when 6.0 M NaOH was added to this same tube and explain.

Part E


Give the results for any Critical Thinking Applications completed and answer all questions given in the Procedure for these exercises.

QUESTIONS

- For each of the following pairs of solutes and solvent, predict whether the solute would be soluble or insoluble. After making your predictions, you can check your answers by looking up the compounds in *The Merck Index* or the *CRC Handbook of Chemistry and Physics*. Generally, *The Merck Index* is the easier reference book to use. If the substance has a solubility greater than 40 mg/mL, you may conclude that it is soluble.
 - Malic acid in water

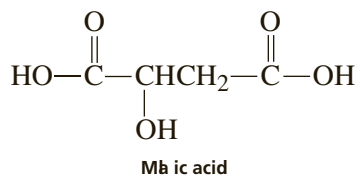


b. Naphthalene in water

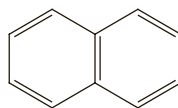


Naphthalene

c. Amphetamine in ethyl alcohol

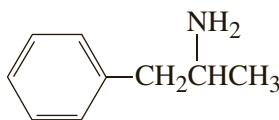


Amphetamine

d. Aspirin in water

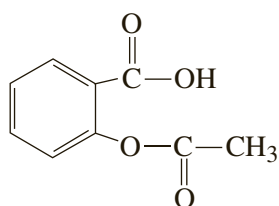


Aspirin

e. Succinic acid in hexane (Note: The polarity of hexane is similar to that of petroleum ether.)

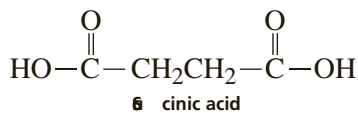


succinic acid

f. Ibuprofen in diethyl ether

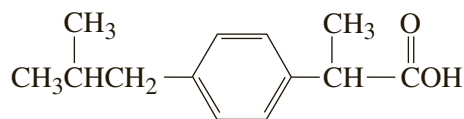


Ibuprofen

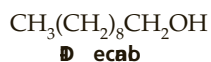g. 1-Decanol (*n*-decyl alcohol) in water


1-decanol

2. Predict whether the following pairs of liquids would be miscible or immiscible:

- Water and methyl alcohol
- Hexane and benzene
- Methylene chloride and benzene
- Water and toluene
- Cyclohexanone and water

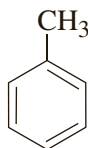


Toluene

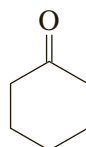


Cyclohexanone


f. Ethyl alcohol and isopropyl alcohol

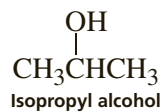

3. Would you expect ibuprofen (see 1f) to be soluble or insoluble in 1.0 M NaOH? Explain.
4. Thymol is very slightly soluble in water and very soluble in 1.0 M NaOH. Explain.

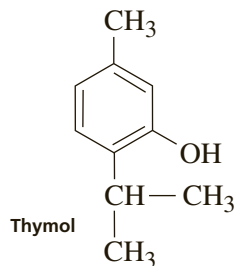

5. Although cannabinol and methyl alcohol are both alcohols, cannabinol is very slightly soluble in methyl alcohol at room temperature. Explain.

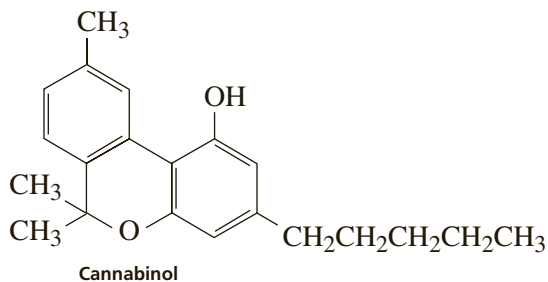

Crystallization

Crystallization

Vacuum filtration

Melting point

Finding a crystallization solvent

Mixture melting point

Critical thinking application

The purpose of this experiment is to introduce the technique of crystallization, the most common procedure used to purify crude solids in the organic laboratory. For a thorough discussion of crystallization, read Technique 11 before proceeding because an understanding of this material is assumed in this experiment.


In Experiments 3A and 3B, you will carry out a crystallization of impure sulfanilamide using 95% ethyl alcohol as the solvent. Sulfanilamide is one of the sulfa drugs, the first generation of antibiotics to be used in successfully treating many major diseases, such as malaria, tuberculosis, and leprosy (see the essay “Sulfa Drugs,” which precedes Experiment 51).

In Experiments 3A and 3B, and in most of the experiments in this textbook, you are told what solvent to use for the crystallization procedure. Some of the factors involved in selecting a crystallization solvent for sulfanilamide are discussed in Technique 11, Section 11.5. The most important consideration is the shape of the solubility curve for the solubility vs. temperature data. As can be seen in Technique 11, Figure 11.2, the solubility curve for sulfanilamide in 95% ethyl alcohol indicates that ethyl alcohol is an ideal solvent for crystallizing sulfanilamide.

The purity of the final material after crystallization will be determined by performing a melting point on your sample. You will also weigh your sample and calculate the percentage recovery. It is impossible to obtain a 100% recovery. This is true for several reasons: There will be some experimental loss, the original sample is not 100% sulfanilamide, and some sulfanilamide is soluble in the solvent even at 0°C. Because of this last fact, some sulfanilamide will remain dissolved in the **mother liquor** (the liquid remaining after crystallization has taken place). Sometimes it is worth isolating a second crop of crystals from the mother liquor, especially if you have performed a synthesis requiring many hours of work and the amount of product is relatively small. This can be accomplished by heating the mother liquor to evaporate some of the solvent and then cooling the resultant solution to induce a second crystallization. The purity of the second crop will not be as good as that of the first crop, however, because the concentration of the impurities will be greater in the mother liquor after some of the solvent has been evaporated.

Two procedures are given here for crystallizing sulfanilamide: a semimicroscale procedure using an Erlenmeyer flask and a Hirsch funnel (Experiment 3A) and a microscale procedure with a Craig tube (Experiment 3B). Your instructor may assign both or just one of these procedures.

In Experiment 3C you will be given an impure sample of the organic compound fluorene (see structure that follows). You will use an experimental procedure for determining which one of three possible solvents is the most appropriate. The three solvents will illustrate three very different solubility behaviors: One of the solvents will be an appropriate solvent for crystallizing fluorene. In a second solvent, fluorene will be highly soluble, even at room temperature. Fluorene will be relatively insoluble in the third solvent, even at the boiling point of the solvent. Your task will be to find the appropriate solvent for crystallization and then perform a crystallization on this sample.

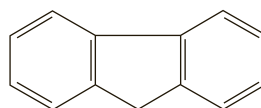

Fluorene

You should be aware that not all crystallizations will look the same. Crystals have many different shapes and sizes, and the amount of mother liquor visible at the end of the crystallization may vary significantly. The crystallizations of sulfanilamide and fluorene will appear significantly different even though the purity of the crystals in each case should be very good.

In Experiment 3D of this experiment, you will determine the identity of an unknown using the melting point technique. The **mixture melting point** technique is introduced in this part.

REQUIRED READING

<i>Review:</i>	Technique 10	Solubility
<i>New:</i>	Technique 8	Filtration, Sections 8.3 and 8.5
	Technique 9	Physical Constants of Solids: The Melting Point
	Technique 11	Crystallization: Purification of Solids

SUGGESTED WASTE DISPOSAL

Dispose of all organic wastes into the nonhalogenated organic waste container.

3A EXPERIMENT 3A

Semimicroscale Crystallization— Erlenmeyer Flask and Hirsch Funnel

This experiment assumes a familiarity with the general semimicroscale crystallization procedure (Technique 11, Section 11.3). In this experiment, Step 2 in Technique 11, Figure 11.3 (removal of insoluble impurities) will not be required. Although the impure sample may have a slight color, it will also not be necessary to use a decolorizing agent (Technique 11, Section 11.7). Leaving out these steps makes the crystallization easier to perform. Furthermore, very few experiments in

this textbook require either of these techniques. If a filtration or decolorizing step is ever required, Technique 11 describes these procedures in detail.

Pre-Lab Calculations

1. Calculate how much 95% ethyl alcohol will be required to dissolve 0.3 g of sulfanilamide at 78°C. Use the data for the graph in Technique 11, Figure 11.2 to make this calculation. The reason for making this calculation is so that you will know ahead of time the approximate amount of hot solvent you will be adding.
2. Using the volume of solvent calculated in Step 1, calculate how much sulfanilamide will remain dissolved in the mother liquor after the mixture is cooled to 0°C.

To dissolve the sulfanilamide in the minimum of hot (boiling or almost boiling) solvent, you must keep the mixture at (or near) the boiling point of 95% ethyl alcohol during the entire procedure. You will likely add more solvent than the amount you calculated because some solvent will evaporate. The amount of solvent is calculated only to indicate the approximate amount of solvent required. You should follow the procedure to determine the correct amount of solvent needed.

PROCEDURE

Preparations

Weigh 0.30 g of impure sulfanilamide¹ and transfer this solid to a 10-mL Erlenmeyer flask. Note the color of the impure sulfanilamide. To a second Erlenmeyer flask, add about 6 mL of 95% ethyl alcohol and a boiling stone. Heat the solvent on a *warm* hot plate until it is boiling.² Because 95% ethyl alcohol boils at a relatively low temperature (78°C), it evaporates rapidly. Setting the temperature of the hot plate too high will result in too much loss of solvent through evaporation.

Dissolving the Sulfanilamide

Before heating the flask containing the sulfanilamide, add enough hot solvent with a Pasteur pipette to barely cover the crystals. Then heat the flask containing the sulfanilamide until the solvent is boiling. At first, this may be difficult to see because so little solvent is present. Add another small portion of solvent (several drops, or about 0.25 mL), continue to heat the flask, and swirl the flask frequently. You may swirl the flask while it is on the hot plate or, for more vigorous swirling, remove it from the hot plate for a few seconds while you swirl it. When you have swirled the flask for 10–15 seconds, check to see if the solid has dissolved. If it has not, add another portion of solvent. Heat the flask again with swirling until the solvent boils. Then swirl the flask for 10–15 seconds, frequently returning the flask to the hot plate so that the temperature of the mixture does not drop. Continue repeating the process of adding solvent, heating, and swirling until all the solid has dissolved completely. Note that it is essential to add just enough solvent to dissolve the solid—neither too much nor too little. Because 95% ethyl alcohol is very volatile, you need to perform this entire procedure fairly rapidly. Otherwise, you may lose solvent nearly as quickly as you are adding it, and this procedure will take a very long time. The time from the first addition of solvent until the solid dissolves completely should not be longer than 10–15 minutes.

¹ The impure sulfanilamide contains 5% fluorenone, a yellow compound, as the impurity.

² To prevent bumping in the boiling solvent, you may want to place a Pasteur pipette in the flask. Use a 25-mL flask so that the Pasteur pipette does not tip the flask over. This is a convenient method because a Pasteur pipette will also be used to transfer the solvent.

Crystallization

Remove the flask from the heat and allow the solution to cool *slowly* (see Technique 11, Section 11.3, Part C, for suggestions). Cover the flask with a small watch glass or stopper the flask. Crystallization should begin by the time the flask has cooled to room temperature. If it has not, scratch the inside surface of the flask with a glass rod (not fire-polished) to induce crystallization (Technique 11, Section 11.8, Part A). When it appears that no further crystallization is occurring at room temperature, place the flask in an ice-water bath using a beaker (Technique 11, Section 11.8). Be sure that both water and ice are present and that the beaker is small enough to prevent the flask from tipping over.

Filtration

When crystallization is complete, vacuum-filter the crystals using a Hirsch funnel (see Technique 8, Section 8.3, and Figure 8.5). (If you will be performing the Optional Exercise at the end of this procedure, you must save the mother liquor from this filtration procedure. Therefore, the filter flask should be clean and dry.) Moisten the filter paper with a few drops of 95% ethyl alcohol and turn on the vacuum (or aspirator) to the fullest extent. Use a spatula to dislodge the crystals from the bottom of the flask before transferring the material to the Hirsch funnel. Swirl the mixture in the flask and pour the mixture into the funnel, attempting to transfer both crystals and solvent. You need to perform these two steps (“swirl and dump”) quickly, before the crystals have completely resettled on the bottom of the flask. (You may need to do this in portions, depending on the size of your Hirsch funnel.) When the liquid has passed through the filter, repeat this procedure until you have transferred all the liquid to the Hirsch funnel. At this point, there will usually be some crystals remaining in the flask. Using your spatula, scrape out as many of the crystals as possible from the flask. Add about 1 mL of *ice-cold* 95% ethyl alcohol (measured with a calibrated Pasteur pipette) to the flask. Swirl the liquid in the flask and then pour the remaining crystals and alcohol into the Hirsch funnel. Not only does this additional solvent help transfer the remaining crystals to the funnel but the alcohol also rinses the crystals already on the funnel. This washing step should be done whether or not it is necessary to use the wash solvent for transferring crystals. If necessary, repeat with another 1-mL portion of ice-cold alcohol. Wash the crystals with a total of about 2 mL of ice-cold solvent.

Continue drawing air through the crystals on the Hirsch funnel by suction for about five minutes. Transfer the crystals onto a preweighed watch glass for air-drying. (Save the mother liquor in the filter flask if you will be doing the Optional Exercise.) Separate the crystals as much as possible with a spatula. The crystals should be completely dried within 10–15 minutes. You can usually determine if the crystals are still wet by observing whether or not they stick to a spatula or stay together in a clump. Weigh the dry crystals and calculate the percent recovery. Compare the color of the pure sulfanilamide to the impure sulfanilamide at the beginning of the experiment. Determine the melting point of the pure sulfanilamide and the original impure material. The literature melting point for pure sulfanilamide is 163–164°C. At the option of the instructor, turn in your crystallized material in a properly labeled container.

Comments on the Crystallization Procedure

1. Do not heat the crude sulfanilamide until you have added some solvent. Otherwise, the solid may melt and possibly form an oil, which may not crystallize easily.
2. When you are dissolving the solid in hot solvent, the solvent should be added in small portions with swirling and heating. The procedure calls for a specific amount (about 0.25 mL), which is appropriate for this experiment. However,

the actual amount you should add each time you perform a crystallization may vary, depending on the size of your sample and the nature of the solid and solvent. You will need to make this judgment when you perform this step.

3. One of the most common mistakes is to add too much solvent. This can happen most easily if the solvent is not hot enough or if the mixture is not stirred sufficiently. If too much solvent is added, the percent recovery will be reduced; it is even possible that no crystals will form when the solution is cooled. If too much solvent is added, you must evaporate the excess by heating the mixture. Using a nitrogen or air stream directed into the container will accelerate the evaporation process (see Technique 7, Section 7.10).
4. Sulfanilamide should crystallize as large, beautiful needles. However, this will not always happen. If the crystals form too rapidly or if there is not enough solvent, they will tend to be smaller, perhaps even appearing as a powder. Furthermore, many substances crystallize in other characteristic shapes, such as plates or prisms.
5. When the solvent is water or when the crystals form as a powder, it will be necessary to dry the crystals longer than 10–15 minutes. Overnight drying may be necessary, especially with water.

Optional Exercise

Transfer the mother liquor to a tared (preweighed) test tube. Place the test tube in a hot water bath and evaporate all the solvent from the mother liquor. Use a stream of nitrogen or air directed into the test tube to speed up the rate of evaporation (see Technique 7, Section 7.10). Cool the test tube to room temperature and dry the outside. Weigh the test tube with solid. Compare this to the weight calculated in the Pre-Lab Calculations. Determine the melting point of this solid and compare it to the melting point of the crystals obtained by crystallization.

3B

EXPERIMENT 3B

Microscale Crystallization—Craig Tube

This experiment assumes familiarity with the general microscale crystallization procedure (Technique 11, Section 11.4). In this experiment, Step 2 in Technique 11, Figure 11.6 (removal of insoluble impurities) will not be required. Although the impure sample may have a slight color, it also will not be necessary to use a decolorizing agent (Technique 11, Section 11.7). Leaving out these steps makes the crystallization easier to perform. Furthermore, very few experiments in this textbook require either of these techniques. When a filtration or decolorizing step is required, Technique 11 describes these procedures in detail.

Pre-Lab Calculations

1. Calculate how much 95% ethyl alcohol will be required to dissolve 0.1 g of sulfanilamide at 78°C. Use the data for the graph in Technique 11, Figure 11.2 to make this calculation. Make this calculation so that you will know the approximate amount of hot solvent you will be adding.
2. Using the volume of solvent calculated in Step 1, calculate how much sulfanilamide will remain dissolved in the mother liquor after the mixture is cooled to 0°C.

To dissolve the sulfanilamide in the minimum of hot (boiling or almost boiling) solvent, you must keep the mixture at (or near) the boiling point of 95% ethyl alcohol during the entire procedure. You will likely add more solvent than the amount you calculated because some solvent will evaporate. Use this calculated amount only as a guide: you should follow the procedure to determine the correct amount of solvent needed.

PROCEDURE

Preparations

Weigh 0.10 g of impure sulfanilamide³ and transfer this solid to a Craig tube. Note the color of the impure sulfanilamide. To a small test tube, add 2–3 mL of 95% ethyl alcohol and a boiling stone. Heat the solvent on a *warm (not hot)* hot plate with an aluminum block until the solvent is boiling.⁴ Setting the temperature of the hot plate too high will result in too much loss of solvent through evaporation.

CAUTION



In performing the following procedure, keep in mind that the mixture in the Craig tube may erupt out of the tube if it becomes superheated. You can prevent this by stirring the mixture constantly with the spatula and by avoiding overheating the mixture.

Dissolving the Sulfanilamide

Before heating the Craig tube containing the sulfanilamide, add enough hot solvent with a Pasteur pipette to barely cover the crystals. Then heat the Craig tube containing the sulfanilamide until the solvent is boiling. At first, this may be difficult to see because so little solvent is present. Add another small portion of solvent (one or two drops), continue to heat the Craig tube, and stir the mixture by rapidly twirling a microspatula between your fingers. When you have stirred the mixture for 10–15 seconds, check to see whether the solid has dissolved. If it has not, add another portion (one or two drops) of solvent. Heat the Craig tube again with stirring until the solvent boils. Then stir the tube for 10–15 seconds. Continue repeating this process of adding solvent, heating, and stirring until all the solid has dissolved completely. Note that it is essential to add just enough solvent to dissolve the solid—neither too much nor too little. Because 95% ethyl alcohol is very volatile, you need to perform this entire procedure fairly rapidly. Otherwise, you may lose solvent nearly as rapidly as you are adding it, and this procedure will take a very long time. The time from the first addition of solvent until the solid dissolves completely should be no longer than 10–15 minutes.

Crystallization

Remove the Craig tube from the heat and insert the inner plug into the opening. Allow the Craig tube to cool slowly to room temperature by placing it into a 10-mL Erlenmeyer flask (see Technique 11, Section 11.4C). Crystallization should begin by the time the Craig tube has cooled to room temperature. If it has not, *gently* scratch the inside surface of the tube with a glass rod (not fire-polished) to induce crystallization

³ See footnote 1 in Experiment 3A.

⁴ You may also use a hot water bath to heat the solvent in the test tube and to heat the Craig tube. The temperature of the water bath should be about 80°C.

(Technique 11, Section 11.7, Part B).⁵ When it appears that no further crystallization is occurring at room temperature, place the Craig tube in an ice-water bath using a beaker (Technique 6, Section 6.5). Be sure that both water and ice are present and that the beaker is small enough to prevent the Craig tube from tipping over.

Isolation of Crystals

When crystallization is complete, place the Craig tube in a centrifuge tube and separate the crystals from the mother liquor by centrifugation. Follow the procedure in Technique 11, Section 11.7.

Using the copper wire, pull the Craig tube out of the centrifuge tube. If the crystals collected on the end of the inner plug, remove the plug and scrape the crystals with a spatula onto a preweighed watch glass for air-drying. Otherwise, it will be necessary to scrape the crystals from the inside surface of the outer part of the Craig tube. If you will be doing the Optional Exercise, save the mother liquor in the centrifuge tube. Separate the crystals as much as possible with a spatula. The crystals should be completely dried within 5–10 minutes. You can determine if the crystals are still wet by observing whether or not they stick to a spatula or stay together in a clump. Weigh the dry crystals and calculate the percent recovery. Compare the color of the pure sulfanilamide to that of the impure sulfanilamide at the beginning of the experiment. Determine the melting point of both the pure sulfanilamide and the original impure material. The literature melting point for pure sulfanilamide is 163–164°C. At the option of the instructor, turn in your crystallized material in a properly labeled container.

For additional information about crystallization, see “Comments on the Crystallization Procedure” in Experiment 3A.

Optional Exercise

See “Optional Exercise” in Experiment 3A.

3C

EXPERIMENT 3C

Selecting a Solvent to Crystallize a Substance

In this experiment you will be given an impure sample of fluorene.⁶ Your goal will be to find a good solvent for crystallizing the sample. You should try water, methyl alcohol, and toluene. After you have determined which is the best solvent, crystallize the remaining material. Finally, determine the melting point of the purified compound and of the impure sample.

PROCEDURE

Selecting a Solvent

Perform the procedure given in Technique 11, Section 11.6 with three separate samples of impure fluorene. Use the following solvents: methyl alcohol, water, and toluene.

⁵ An alternative method for inducing crystallization is to dip a microspatula into the solution. Then allow the solvent to evaporate so that a small amount of solid will form on the surface of the spatula. When placed back into the solution, the solid will seed the solution.

⁶ The impure fluorene contains 5% fluorenone, a yellow compound.

Crystallizing the Sample

After you have found a good solvent, crystallize the impure fluorene using a semimicroscale (Erlenmeyer flask and Hirsch funnel) or a microscale (Craig tube) procedure. Use 0.3 g of impure fluorene if you follow the semimicroscale procedure, or use 0.05 g if you follow the microscale procedure. Weigh the impure sample carefully, and be sure to keep a little of the impure sample on which to perform a melting point. If you perform a semimicroscale crystallization, you may need to use a size of Erlenmeyer flask different from the one specified in the procedure. This decision should be made based on the amount of sample you will be crystallizing and how much solvent you think will be needed. Transfer the crystals to a preweighed watch glass and allow them to air-dry. If water was used as the solvent, you may need to let the crystals sit out overnight for drying because water is less volatile than most organic solvents. Weigh the dried sample and calculate the percent recovery. Determine the melting point of both the pure sample and the original impure material. The literature melting point for pure fluorene is 116–117°C. At the option of the instructor, turn in your crystallized material in a properly labeled container.

3D EXPERIMENT 3D

Mixture Melting Points

In Experiments 3A and 3B of this experiment, the melting point was used to determine the purity of a known substance. In some situations the melting point can also be used to determine the identity of an unknown substance.

In Experiment 3D, you will be given a pure sample of an unknown from the following list:

Compound	Melting Point (°C)
Acetylsalicylic acid	138–140
Benzoic acid	121–122
Benzoin	135–136
Dibenzoyl ethylene	108–111
Succinimide	122–124
<i>o</i> -Toluic acid	108–110

Your goal is to determine the identity of the unknown using the melting-point technique. If all of the compounds in the list had distinctly different melting points, it would be possible to determine the identity of the unknown by just taking its melting point. However, each of the compounds in this list has a melting point that is close to the melting point of another compound in the list. Therefore, the melting point of the unknown will allow you to narrow down the choices to two compounds. To determine the identity of your compound, you must perform mixture melting points of your unknown and each of the two compounds with similar melting points. A mixture melting point that is depressed and has a wide range indicates that the two compounds in the mixture are different.

PROCEDURE

Obtain an unknown sample and determine its melting point. Determine mixture melting points (see Technique 9, Section 9.4) of your unknown and all compounds from the previous list that have similar melting points. To prepare a sample for a mixture melting point, use a spatula or a glass stirring rod to grind equal amounts of your unknown and the known compound in a watch glass. Record all melting points and state the identity of your unknown.


3E

EXPERIMENT 3E

Critical Thinking Application

The goal of the exercise is to find an appropriate solvent to crystallize a given compound. Rather than doing this experimentally, you will try to predict which one of three given solvents is the best. For each compound, one of the solvents has the desired solubility characteristics to be a good solvent for crystallization. In a second solvent, the compound will be highly soluble, even at room temperature. The compound will be relatively insoluble in the third solvent, even at the boiling point of the solvent. After making your predictions, you will check them by looking up the appropriate information in *The Merck Index*.

For example, consider naphthalene, which has the following structure:

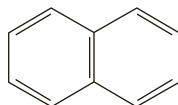

Naphthalene

Consider the three solvents ether, water, and toluene. (Look up their structures if you are unsure. Remember that ether is also called diethyl ether.) Based on your knowledge of polarity and solubility behavior, make your predictions. It should be clear that naphthalene is insoluble in water because naphthalene is a hydrocarbon that is nonpolar and water is very polar. Both toluene and ether are relatively nonpolar, so naphthalene is most likely soluble in both of them. One would expect naphthalene to be more soluble in toluene because both naphthalene and toluene are hydrocarbons. In addition, they both contain benzene rings, which means that their structures are very similar. Therefore, according to the solubility rule “Like dissolves like,” one would predict that naphthalene is very soluble in toluene. Perhaps it is too soluble in toluene to be a good crystallizing solvent. If so, then ether would be the best solvent for crystallizing naphthalene.

These predictions can be checked with information from *The Merck Index*. Finding the appropriate information can be somewhat difficult, especially for beginning organic chemistry students. Look up *naphthalene* in *The Merck Index*. The entry for *naphthalene* states, “Monoclinic prismatic plates from ether.”

This statement means that naphthalene can be crystallized from ether. It also gives the type of crystal structure. Unfortunately, sometimes the crystal structure is given without reference to the solvent. Another way to determine the best solvent is by looking at solubility-vs.-temperature data. A good solvent is one in which the solubility of the compound increases significantly as the temperature increases. To determine whether the solid is too soluble in the solvent, check the solubility at room temperature. In Technique 11, Section 11.6, you were instructed to add 0.5 mL of solvent to 0.05 g of compound. If the solid completely dissolved, then the solubility at room temperature was too great. Follow this same guideline here. For naphthalene, the solubility in toluene is given as 1 g in 3.5 mL. When no temperature is given, room temperature is understood. By comparing this to the 0.05 g in 0.5 mL ratio, it is clear that naphthalene is too soluble in toluene at room temperature for toluene to be a good crystallizing solvent. Finally, *The Merck Index* states that naphthalene is insoluble in water. Sometimes no information is given about solvents in which the compound is insoluble. In that case, you would rely on your understanding of solubility behavior to confirm your predictions.

When using *The Merck Index*, you should be aware that alcohol is listed frequently as a solvent. This generally refers to 95% or 100% ethyl alcohol. Because 100% (absolute) ethyl alcohol is more expensive than 95% ethyl alcohol, the cheaper grade is usually used in the chemistry lab. Finally, benzene is frequently listed as a solvent. Because benzene is a known carcinogen, it is rarely used in student labs. Toluene is a suitable substitute; the solubility behavior of a substance in benzene and toluene is so similar that you may assume any statement made about benzene also applies to toluene.


Exercise. For each of the following sets of compounds (the solid is listed first, followed by the three solvents), use your understanding of polarity and solubility to predict

1. The best solvent for crystallization
2. The solvent in which the compound is too soluble
3. The solvent in which the compound is not sufficiently soluble

Then check your predictions by looking up each compound in *The Merck Index*.


1. Phenanthrene; toluene, 95% ethyl alcohol, water

Phenanthrene

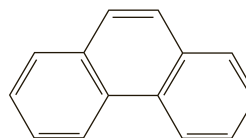


2. Cholesterol; ether, 95% ethyl alcohol, water


Cholesterol

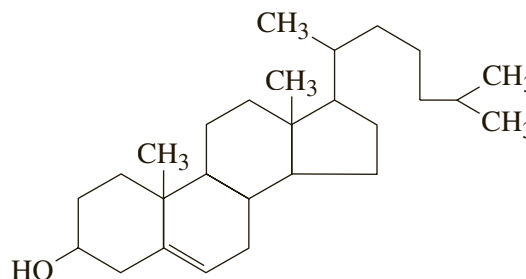

3. Acetaminophen; toluene, 95% ethyl alcohol, water

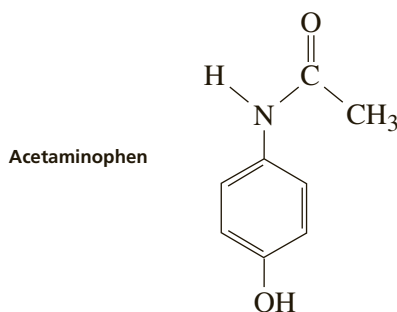

4. Urea; hexane, 95% ethyl alcohol, water

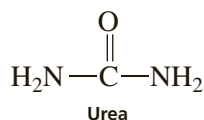

REPORT
Experiments 3A and 3B

1. Report the melting points for both the impure sulfanilamide and the crystallized sulfanilamide and comment on the differences. Also, compare these to the literature value. Based on the melting point of the crystallized sulfanilamide, is it pure? Also comment on the purity based on the color of the crystals. Report both the original weight of the impure sulfanilamide and the weight of the crystallized sulfanilamide. Calculate the percentage recovery and comment on several sources of loss.
2. If you completed the Optional Exercise (isolating the solid dissolved in the mother liquor), do the following:
 - a. Make a table with the following information:
 - i. Weight of impure sulfanilamide used in the crystallization procedure
 - ii. Weight of pure sulfanilamide after crystallization
 - iii. Weight of sulfanilamide plus impurity recovered from the mother liquor (see Experiment 3A or 3B, Optional Exercise)
 - iv. Total of items ii and iii (total weight of sulfanilamide plus impurity isolated)
 - v. Calculated weight of sulfanilamide in the mother liquor (see Experiment 3A or 3B, Pre-Lab Calculations)
 - b. Comment on any differences between the values in items i and iv. Should they be the same? Explain.
 - c. Comment on any differences between items iii and v. Should they be the same? Explain.
 - d. Report the melting point of the solid recovered from the mother liquor. Compare this to the melting points of the crystallized sulfanilamide. Should they be the same? Explain.

Experiment 3C

1. For each of the three solvents (methyl alcohol, water, and toluene), describe the results from the tests for selecting a good crystallizing solvent for fluorene.

Explain these results in terms of polarity and solubility predictions (see “Guidelines for Predicting Polarity and Solubility,” in Technique 10, Section 10.2A).

- Report the melting points for both the impure fluorene and the crystallized fluorene and comment on the differences. What is the literature value for the melting point of fluorene? Report the original weight of both the impure fluorene and the weight of the crystallized fluorene. Calculate the percentage recovery and comment on several sources of loss.
- The solubility of fluorene in each solvent used in Experiment 3B corresponds to one of the three curves shown in Technique 11, Figure 11.1. For each solvent, indicate which curve best describes the solubility of fluorene in that solvent.

Experiment 3D

Record all melting points and state the identity of your unknown.

Experiment 3E

For each compound assigned, state your predictions, along with an explanation. Then give the relevant information from *The Merck Index* that supports or contradicts your predictions. Try to explain any differences between your predictions and information found in *The Merck Index*.

QUESTIONS

- Consider a crystallization of sulfanilamide in which 10 mL of hot 95% ethyl alcohol is added to 0.10 g of impure sulfanilamide. After the solid has dissolved, the solution is cooled to room temperature and then placed in an ice-water bath. No crystals form, even after scratching with a glass rod. Explain why this crystallization failed. What would you have to do at this point to make the crystallization work? Assume that starting over again with a new sample is not an option. (You may need to refer to Technique 11, Figure 11.2.)
- Benzyl alcohol (bp 205°C) was selected by a student to crystallize fluorenol (mp 153–154°C) because the solubility characteristics of this solvent are appropriate. However, this solvent is not a good choice. Explain.
- A student performs a crystallization on an impure sample of biphenyl. The sample weighs 0.5 g and contains about 5% impurity. Based on his knowledge of solubility, the student decides to use benzene as the solvent. After crystallization, the crystals are dried and the final weight is found to be 0.02 g. Assume that all steps in the crystallization are performed correctly, there are no spills, and the student lost very little solid on any glassware or in any of the transfers. Why is the recovery so low?

Extraction

Extraction

Critical thinking application

Extraction is one of the most important techniques for isolating and purifying organic substances. In this method, a solution is mixed thoroughly with a second solvent that is **immiscible** with the first solvent. (Remember that immiscible liquids do not mix; they form two phases, or layers.) The solute is extracted from one solvent into the other because it is more soluble in the second solvent than in the first.

The theory of extraction is described in detail in Technique 12, Sections 12.1–12.2. You should read these sections before continuing this experiment. Because solubility is the underlying principle of extraction, you may also wish to reread the introduction to the experiment on solubility.

Extraction is not only a technique used by organic chemists but it is also used to produce common products with which you are familiar. For example, vanilla extract, the popular flavoring agent, was originally extracted from vanilla beans using alcohol as the organic solvent. Decaffeinated coffee is made from coffee beans that have been decaffeinated by an extraction technique (see the essay, “Caffeine,” which precedes Experiment 13). This process is similar to the procedure in Experiment 4A of this experiment, in which you will extract caffeine from an aqueous solution.

The purpose of this experiment is to introduce the microscale technique for performing extractions and allow you to practice this technique. This experiment also demonstrates how extraction is used in organic experiments.

REQUIRED READING

<i>New:</i>	Technique 12	Extraction
	Essay	Caffeine
<i>Review:</i>	Technique 10	Solubility

SPECIAL INSTRUCTIONS

Be careful when handling methylene chloride. It is a toxic solvent, and you should not breathe its fumes excessively or spill it on yourself.

In Experiment 4B, it is advisable to pool the data for the distribution coefficients and calculate class averages. This will compensate for differences in the values due to experimental error.

SUGGESTED WASTE DISPOSAL

You must dispose of all methylene chloride in a waste container marked for the disposal of halogenated organic wastes. Place all other organic wastes into the nonhalogenated organic waste container. The aqueous solutions obtained after the extraction steps must be disposed of in the container designated for aqueous waste.

4A EXPERIMENT 4A

Extraction of Caffeine

One of the most common extraction procedures involves using an organic solvent (nonpolar or slightly polar) to extract an organic compound from an aqueous solution. Because water is highly polar, the mixture will separate into two layers, or phases: an aqueous layer and an organic (nonpolar) layer.

In this experiment, you will extract caffeine from an aqueous solution using methylene chloride. You will perform the extraction step three times using three separate portions of methylene chloride. Because methylene chloride is more dense than water, the organic layer (methylene chloride) will be on the bottom. After each extraction, you will remove the organic layer. The organic layers from all three extractions will be combined and dried over anhydrous sodium sulfate. After transferring the dried solution to a preweighed container, you will evaporate the methylene chloride and determine the weight of caffeine extracted from the aqueous solution. This extraction procedure succeeds because caffeine is much more soluble in methylene chloride than in water.

Pre-Lab Calculation

In this experiment, 0.070 g of caffeine is dissolved in 4.0 mL of water. The caffeine is extracted from the aqueous solution three times with 2.0-mL portions of methylene chloride. Calculate the total amount of caffeine that can be extracted into the three portions of methylene chloride (see Technique 12, Section 12.2). Caffeine has a distribution coefficient of 4.6, between methylene chloride and water.

PROCEDURE

NOTE: To obtain good results, you should make all weighings as accurately as possible, preferably on a balance that is accurate to within 0.001 g.

Preparation

Before beginning this experiment, check your screw-cap centrifuge tube for leaks.¹ Add exactly 0.070 g of caffeine to the centrifuge tube. Then add 4.0 mL of water to the tube. Cap the tube and shake it vigorously for several minutes until the caffeine dissolves completely. It may be necessary to heat the mixture slightly to dissolve all the caffeine.

¹ Place about 2 mL of water in the tube. Cap it and shake vigorously. If it leaks, try screwing the cap on more tightly or use a different cap. Sometimes you may need to replace the centrifuge tube itself. Discard the water in the tube.

Extraction

Add 2.0 mL of methylene chloride to the tube. The two layers must be mixed thoroughly so that as much caffeine as possible is transferred from the aqueous layer to the methylene chloride layer. However, if the mixture is mixed too vigorously, it may form an emulsion. Emulsions look like a third frothy layer between the other two layers, and they can make it difficult for the layers to separate. The best way to prevent an emulsion is to shake gently at first and observe whether the layers separate. If they separate quickly, continue to shake, but now more vigorously. The correct way to shake is to invert the tube and right it in a rocking motion. A good rate of shaking is about one rock per second. When it is clear that an emulsion is not forming, you may shake it more vigorously, perhaps two to three times per second. (Note that it is usually not prudent to shake the heck out of it!) Shake the tube for about one minute.

After shaking, place the tube in a test tube rack or beaker and let it stand until the layers separate completely.² It may be necessary to tap the sides of the tube to force all the methylene chloride layer to the bottom of the vial. Occasionally, a drop of water will get stuck in the very bottom part of the tube, below the methylene chloride layer. If this happens, depress the bulb slightly and try to draw the water drop into a Pasteur pipette. Transfer this drop to the upper layer.

Using a Pasteur pipette, you should now transfer the organic (bottom) layer into a test tube. Ideally, the goal is to remove all the organic layer without transferring any of the aqueous layer. However, this is difficult to do. Try to squeeze the bulb so that when it is released completely, you will draw up the amount of liquid that you desire. If you have to hold the bulb in a partially depressed position while making a transfer, it is likely that you will spill some liquid. It is also best to transfer the liquid in two steps. First, depress the bulb so that most (about 75%) of the bottom layer will be drawn into the pipette. Place the tip of the pipette squarely in the V at the bottom of the centrifuge tube and release the bulb slowly. When making the transfer, it is essential that the centrifuge tube and the test tube be held next to each other. A good technique for this is illustrated in Figure 12.6. After transferring the first portion, depress the bulb partially, just enough to draw up the remaining liquid in the bottom layer, and place the tip of the pipette in the bottom of the tube. Draw the liquid into the pipette and transfer this liquid to the test tube.

Repeat this extraction two more times using 2 mL of fresh methylene chloride each time. Combine the organic layer from each of these extractions with the methylene chloride solution from the first extraction.

Drying the Organic Layers

Dry the combined organic layers over granular anhydrous sodium sulfate, following the instructions given in Technique 12, Section 12.9, "Drying Procedure with Anhydrous Sodium Sulfate". Read these instructions carefully and complete Steps 1–3 in the "Microscale Drying Procedure." Step 4 is described in the next section, "Evaporation of Solvent."

Evaporation of Solvent

Transfer the dried methylene chloride solution with a clean, dry Pasteur pipette to a dry, preweighed 10-mL Erlenmeyer flask or test tube while leaving the drying agent behind.³ (If you had to add more than 3–4 microspatulafuls of anhydrous sodium

² If an emulsion has formed, the two layers may not separate on standing. If they do not separate after about 1–2 minutes, it will be necessary to centrifuge the mixture to break the emulsion. Remember to balance the centrifuge by placing a tube of equal weight on the opposite side.

³ It is easier to avoid transferring any drying agent if you use a filter-tip pipette (Technique 8, Section 8.6).


sulfate, rinse the sodium sulfate with about 0.5 mL of fresh methylene chloride. Stir this with a dry spatula and then transfer this solution to the same preweighed flask.) Evaporate the methylene chloride by heating the flask in a hot water bath at about 45°C. This should be done in a hood and can be accomplished more rapidly if a stream of dry air or nitrogen gas is directed at the surface of the liquid (see Technique 7, Section 7.10). When the solvent is evaporated, remove the flask from the bath and dry the outside of the flask. When the flask has cooled to room temperature, weigh it to determine the amount of caffeine that was in the methylene chloride solution. Compare this weight with the amount of caffeine calculated in the Pre-Lab Calculation.

4B EXPERIMENT 4B


Distribution of a Solute between Two Immiscible Solvents

In this experiment, you will investigate how several different organic solids distribute themselves between water and methylene chloride. A solid compound is mixed with the two solvents until equilibrium is reached. The organic layer is removed, dried over anhydrous sodium sulfate, and transferred to a tared container. After evaporating the methylene chloride, the weight of the organic solid that was in the organic layer is determined. By finding the difference, the amount of solute in the aqueous layer can also be determined. The distribution coefficient of the solid between the two layers can then be calculated and related to the polarity of the solid and the polarities of the two liquids.


Three different compounds will be used: benzoic acid, succinic acid, and sodium benzoate. Their structures are given below. You should perform this experiment on one of the solids and share your data with two other students who worked with the other two solids. Alternatively, data from the entire class may be pooled and averaged. Before performing this experiment, it would be helpful if you predict the relative solubilities of the three compounds in the two solvents.

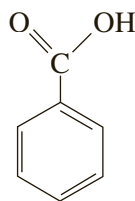

Benzoic acid

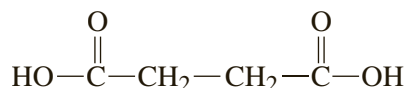

Succinic acid

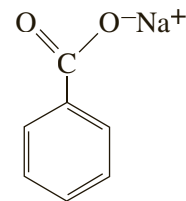

Sodium benzoate

PROCEDURE

NOTE: To obtain good results, you should make all weighings as accurately as possible, preferably on a balance that is accurate to within 0.001g.

Place 0.050 g of one of the solids (benzoic acid, succinic acid, or sodium benzoate) into a 5-mL conical vial. Add 2.0 mL of methylene chloride and 2.0 mL of water to the vial. Cap the vial and shake it as described in Experiment 4A for about 1 minute. Check for undissolved solid. Continue shaking the vial until all the solid is dissolved. After the layers have separated, transfer the bottom organic layer to another vial or a small test tube.

Using the same procedure just described in Experiment 4A (see the section on “Drying the Organic Layers”), dry this organic layer over granular anhydrous sodium sulfate.

Transfer the dried methylene chloride solution with a clean, dry Pasteur pipette to a dry, preweighed test tube, leaving the drying agent behind. Evaporate the methylene chloride by heating the tube in a hot water bath while directing a stream of dry air or nitrogen gas at the surface of the liquid. When the solvent is evaporated, remove the tube from the bath and dry the outside of the tube. When the tube has cooled to room temperature, weigh the tube to determine the amount of solid solute that was in the methylene chloride layer. Determine by difference the amount of the solid that was dissolved in the aqueous layer. Calculate the distribution coefficient for the solid between methylene chloride and water. Because the volume of methylene chloride and water was the same, the distribution coefficient can be calculated by dividing the weight of solute in methylene chloride by the weight of solute in water.

Optional Exercise

Repeat the preceding procedure using 0.050 g of caffeine, 2.0 mL of methylene chloride, and 2.0 mL of water. Determine the distribution coefficient for caffeine between methylene chloride and water. Compare this to the literature value of 4.6.

4C

EXPERIMENT 4C

How Do You Determine Which One Is the Organic Layer?

A common problem that you might encounter during an extraction procedure is not knowing for sure which layer is organic and which is the aqueous one. Although the procedures in this textbook often indicate the expected relative positions of the two layers, not all procedures will give this information, and you should be prepared for surprises. Sometimes knowing the densities of the two solvents is not sufficient, because dissolved substances can significantly increase the density of a solution. It is very important to know the location of the two layers because usually one layer contains the desired product and the other layer is discarded. A mistake at this point in an experiment would be disastrous!

The purpose of this experiment is to give you some practice in determining which layer is aqueous and which layer is organic (see Technique 12, Section 12.8). As described in Section 12.8, one effective technique is to add a few drops of water to each layer after the layers have been separated. If the layer is water, then the drops of added water will dissolve in the aqueous layer and increase its volume. If the added water forms droplets or a new layer, then it is the organic layer.

PROCEDURE

Obtain three test tubes, each containing two layers.⁴ For each tube, you will be told the identity of the two layers, but you will not be told their relative positions. Determine experimentally which layer is organic and which layer is aqueous. Dispose of all these

⁴The three mixtures will likely be (1) water and *n*-butyl chloride, (2) water and *n*-butyl bromide, and (3) *n*-butyl bromide and saturated aqueous sodium bromide.


mixtures into the waste container designated for halogenated organic wastes. After determining the layers experimentally, look up the densities of the various liquids in a handbook to see if there is a correlation between the densities and your results.

4D EXPERIMENT 4D

Use of Extraction to Isolate a Neutral Compound from a Mixture Containing an Acid or Base Impurity

In this experiment you will be given a solid sample containing an unknown neutral compound and an acid or base impurity. The goal is to remove the acid or base by extraction and isolate the neutral compound. By taking the melting point of the neutral compound, you will identify it from a list of possible compounds. There are many organic reactions in which the desired product, a neutral compound, is contaminated by an acid or base impurity. This experiment illustrates how extraction is used to isolate the product in this situation.

In Technique 10, "Solubility," you learned that organic acids and bases can become ions in acid–base reactions (see Section 10.2B, "Solutions in Which the Solute Ionizes and Dissociates"). Before reading on, review this material if necessary. Using this principle, it is possible to separate an acid or base impurity from a neutral compound. The following scheme, which shows how both an acid and a base impurity are removed from the desired product, illustrates how this is accomplished:

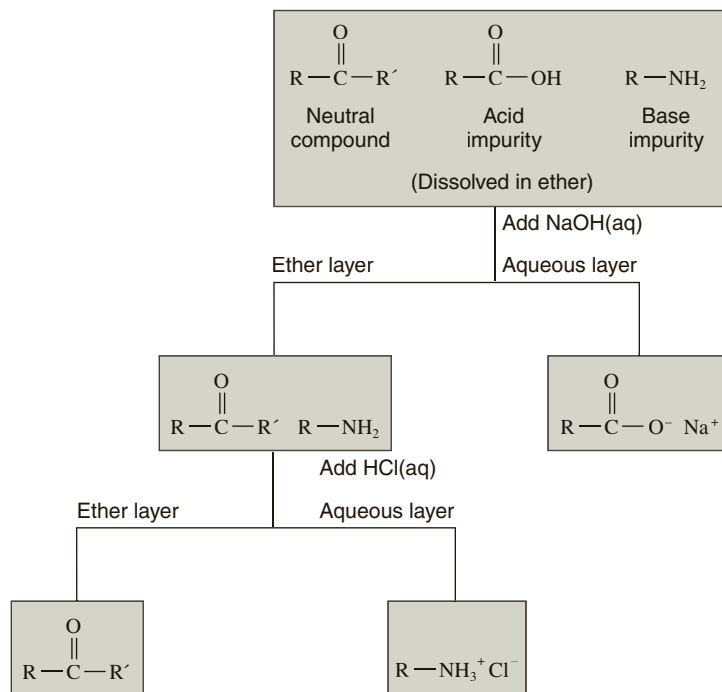

Flow chart showing how acid and base impurities are removed from the desired product.

The neutral compound can now be isolated by removing the water dissolved in the ether and evaporating the ether. Because ether dissolves a relatively large

quantity of water (1.5%), the water must be removed in two steps: In the first step, the ether solution is mixed with a saturated aqueous NaCl solution. Most of the water in the ether layer will be transferred to the aqueous layer in this step (see Technique 12, Section 12.9). Finally, the remainder of the water is removed by drying the ether layer over anhydrous sodium sulfate. The neutral compound can then be isolated by evaporating the ether. In most organic experiments that use a separation scheme such as this, it would be necessary to perform a crystallization step to purify the neutral compound. However, in this experiment the neutral compound should be sufficiently pure at this point to identify it by melting point.

The organic solvent used in this experiment is ether. Recall that the full name for ether is diethyl ether. Because ether is less dense than water, this experiment will give you practice in performing extractions where the nonpolar solvent is less dense than water.

The following procedure details the removal of an acid impurity from a neutral compound and isolating the neutral compound. It contains an additional step that is not normally part of this kind of separation scheme: The aqueous layers from each extraction are segregated and acidified with aqueous HCl. The purpose of this step is to verify that the acid impurity has been removed completely from the ether layer. In the Optional Exercise, the sample contains a neutral compound with a base impurity; however, a detailed procedure is not given. If you are assigned this exercise, you must create a procedure by using the principles discussed in this introduction and by studying the following procedure for isolating the neutral compound from an acid impurity.

PROCEDURE

Isolating a Neutral Compound from a Mixture Containing an Acid Impurity. Add 0.150 g of an unknown mixture⁵ to a screw-cap centrifuge tube. Add 4.0 mL of ether to the tube and cap it. Shake the tube until all the solid dissolves completely.

Add 2.0 mL of 1.0 M NaOH to the tube and shake for 30 seconds. Let the layers separate. Remove the bottom (aqueous) layer, and place this in a test tube labeled "1st NaOH extract." Add another 2.0-mL portion of 1.0 M NaOH to the centrifuge tube and shake for 30 seconds. When the layers have separated, remove the aqueous layer and put this in a test tube labeled "2nd NaOH extract."

With stirring, add 6 M HCl dropwise to each of the two test tubes containing the NaOH extracts until the mixture is acidic. Test the mixture with litmus or pH paper to determine when it is acidic. Observe the amount of precipitate that forms. What is the precipitate? Does the amount of precipitate in each tube indicate that all the acid impurity has been removed from the ether layer containing the unknown neutral compound?

The drying procedure for an ether layer requires the following additional step compared to the procedure for drying a methylene chloride layer (see Technique 12, Section 12.9, "Saturated Salt Solution"). To the ether layer in the centrifuge tube, add 2.0 mL of saturated aqueous sodium chloride. Shake for 30 seconds and let the layers separate. Remove and discard the aqueous layer. With a clean, dry Pasteur pipette, transfer the ether layer (without any water) to a clean, dry test tube. Now dry the ether layer over granular anhydrous sodium sulfate (see Technique 12,

⁵ The mixture contains 0.100 g of one of the neutral compounds given in the list in the following table and 0.050 g of benzoic acid, the acid impurity.

Section 12.9, “Drying Procedure with Anhydrous Sodium Sulfate”). Complete Steps 1–3 in the “Microscale Drying Procedure.” Step 4 is described in the next paragraph.

Transfer the dried ether solution with a clean, dry Pasteur pipette to a dry, preweighed test tube, leaving the drying agent behind. Evaporate the ether by heating the tube in a hot water bath. This should be done in a hood and can be accomplished more rapidly if a stream of dry air or nitrogen gas is directed at the surface of the liquid (see Technique 7, Section 7.10). When the solvent has evaporated, remove the test tube from the bath and dry the outside of the tube. Once the tube has cooled to room temperature, weigh it to determine the amount of solid solute that was in the ether layer. Obtain the melting point of the solid and identify it from the following list:

	Melting Point
Fluorenone	82–85°C
Fluorene	116–117°C
1,2,4,5-Tetrachlorobenzene	139–142°C
Triphenylmethanol	162–164°C

Optional Exercise: Isolating a Neutral Compound from a Mixture Containing a Base Impurity. Obtain 0.150 g of an unknown mixture containing a neutral compound and a base impurity.⁶ Develop a procedure for isolating the neutral compound, using the preceding procedure as a model. After isolating the neutral compound, obtain the melting point and identify it from the list of compounds given above.

4E EXPERIMENT 4E

Critical Thinking Application

PROCEDURE

1. Add 4 mL of water and 2 mL of methylene chloride to a screw-capped centrifuge tube.
2. Add 4 drops of solution A to the centrifuge tube. Solution A is a dilute aqueous solution of sodium hydroxide containing an organic compound.⁷ Shake the mixture for about 30 seconds, using a rapid rocking motion. Describe the color of each layer (see the following table).
3. Add 2 drops of 1 M HCl. Let the solution sit for 1 minute and note the color change. Then shake for about 1 minute, using a rapid rocking motion. Describe the color of each layer.

⁶ The mixture contains 0.100 g of one of the neutral compounds given in the list above and 0.050 g of ethyl 4-aminobenzoate, a base impurity.

⁷ Solution A: Mix 25 mg of 2,6-dichloroindophenol (sodium salt) with 50 mL of water and 1 mL of 1 M NaOH. This solution should be prepared the same day it is used.

- Add 4 drops of 1 M NaOH and shake again for about 1 minute. Describe the color of each layer.

		Color
Step 2	Aqueous	
	Methylene chloride	
Step 3	Aqueous	
	Methylene chloride	
Step 4	Aqueous	
	Methylene chloride	

REPORT

Experiment 4A

- Show your calculations for the amount of caffeine that should be extracted by the three 2.0-mL portions of methylene chloride (see Pre-Lab Calculation).
- Report the amount of caffeine isolated. Compare this weight with the amount of caffeine calculated in the Pre-Lab Calculation. Comment on the similarity or difference.

Experiment 4B

- Report in table form the distribution coefficients for the three solids: benzoic acid, succinic acid, and sodium benzoate.
- Is there a correlation between the values of the distribution coefficients and the polarities of the three compounds? Explain.
- If you completed the Optional Exercise, compare the distribution coefficient you obtained for caffeine with the corresponding literature value. Comment on the similarity or difference.

Experiment 4C

- For each of the three mixtures, report which layer was on the bottom and which one was on the top. Explain how you determined this for each mixture.
- Record the densities for the liquids given in a handbook.
- Was there a correlation between the densities and your results? Explain.

Experiment 4D

- Answer the following questions about the first and second NaOH extracts.
 - Comment on the amount of precipitate for both extracts when HCl is added.
 - What is the precipitate formed when HCl is added?
 - Does the amount of precipitate in each tube indicate that all the acid impurity has been removed from the ether layer containing the unknown neutral compound?
- Report the melting point and weight of the neutral compound you isolated.
- Based on the melting point, what is the identity of this compound?
- Calculate the percent recovery for the neutral compound. List possible sources of loss.

If you completed the Optional Exercise, complete Steps 1–4 for Experiment 4D.

Experiment 4E

Describe fully what occurred in Steps 2, 3, and 4. For each step, include (1) the nature (cation, anion, or neutral species) of the organic compound, (2) an explanation for all the color changes, and (3) an explanation for why each layer is colored

as it is. Your explanation for (3) should be based on solubility principles and the polarities of the two solvents. (*Hint*: It may be helpful to review the sections in your general chemistry textbook that deal with acids, bases, and acid–base indicators.)

REFERENCE

Kelly, T. R. A Simple, Colorful Demonstration of Solubility and Acid/Base Extraction. *Journal of Chemical Education*, 70 (1993): 848.

QUESTIONS

1. Perform an online search using your browser (Mozilla Firefox or Internet Explorer) and a search engine (Google or Bing) to find the structures of the compounds in the questions that follow. Provide a method for separating the mixtures of two or three compounds, dissolved in a solvent diethyl ether. In each case one of the components will be a neutral compound. You should give your answer in the form of a flow chart (see Section 12.12).
 - a. Benzophenone and tributylamine
 - b. 4-Bromoaniline, 3-nitrobenzoic acid, and benzoin
 - c. Fluorenone, octanoic acid, and dicyclohexylamine
 - d. 1-Hexanol and 4-bromoaniline
2. Describe how you could separate and purify compound A from a mixture of two neutral compounds (A and B) when A comprises 95% of the total and B the other 5% of the total. Assume that A and B have similar polarities.
3. Consider a mixture containing 0.5 g of benzil and 0.05 g of benzoin. Your task is to isolate benzil in a pure form. Could you accomplish this using an extraction procedure? If yes, explain how you would do this. If not, suggest another technique that would accomplish this.

A Separation and Purification Scheme

Extraction

Crystallization

Devising a procedure

Critical thinking application

There are many organic experiments in which the components of a mixture must be separated, isolated, and purified. Although detailed procedures are usually given for carrying this out, devising your own scheme can help you understand these techniques more thoroughly. In this experiment, you will devise a separation and purification scheme for a three-component mixture that will be assigned to you. The mixture will contain a neutral organic compound and either an organic acid or base in nearly equal amounts. The third component, also a neutral compound, will be present in a much smaller amount. Your goal will be to isolate in pure form *two* of the three compounds. The components of your mixture may be separated and purified by a combination of acid-base extractions and crystallizations. You will be told the composition of your mixture well in advance of the laboratory period so that you will have time to write a procedure for this experiment.

REQUIRED READING

Review: Technique 11
Technique 12

Crystallization: Purification of Solids
Extractions, Separations, and Drying Agents

SUGGESTED WASTE DISPOSAL

Dispose of all filtrates that may contain 1,4-dibromobenzene or methylene chloride into the container designated for halogenated organic wastes. All other filtrates may be disposed of into the container for nonhalogenated organic wastes.

NOTES TO THE INSTRUCTOR

Students must be told the composition of their mixture well in advance of the laboratory period so that they have enough time to devise a procedure. It is advisable to require that students turn in a copy of their procedure at the beginning of the lab period. You may wish to allow enough time for students to repeat the experiment if their procedure doesn't work the first time or if they want to improve on their percentage recovery and purity. If you allow enough time for students to perform this experiment just once, it will be helpful to put out pure samples of the compounds

Experiment 5 is based on a similar experiment developed by James Patterson, North Seattle Community College, Seattle.

in the mixtures so students can try out different solvents to determine a good solvent for crystallizing each compound.

PROCEDURE

Advance Preparation

Each student will be assigned a mixture of three compounds.¹ Before coming to the laboratory, you must work out a detailed procedure that can be used to separate, isolate, and purify *two* of the compounds in your mixture. You may not be able to specify all the reagents or the volumes required ahead of time, but the procedure should be as complete as possible. It will be helpful to consult the following experiments and techniques:

Experiment 2, "Solubility," Part D

Experiment 4D, "Extraction"

Technique 10, Section 10.2B

Technique 12, Sections 12.11, and 12.12

Keep in mind that the overall purpose of Experiment 4D is somewhat different than your goal in Experiment 5. In Experiment 4D, an acid or base impurity is removed from the neutral compound, but the acid or base impurity is not isolated. In Experiment 5, your separation scheme may include the isolation of an acid or base compound, depending on the composition of your mixture. For this purpose, you will likely find Technique 12, Sections 12.11 and 12.12, more helpful in devising a procedure than Experiment 4D. As part of your advance preparation, you should outline the separation scheme using a flow chart (see Technique 12, Section 12.12).

The following reagents will be available: 1M NaOH, 6M NaOH, 1M HCl, 6M HCl, 1M NaHCO₃, saturated sodium chloride, diethyl ether, 95% ethanol, methanol, isopropyl alcohol, acetone, hexane, toluene, methylene chloride, and granular anhydrous sodium sulfate. Other solvents that can be used for crystallization may also be available.

Separation

The first step in your procedure should be to dissolve about 0.5 g (record exact weight) of the mixture in the minimum amount of diethyl ether or methylene chloride. If more than about 4 mL of a solvent is required, you should use the other solvent. Most of the compounds in the mixtures are more soluble in methylene chloride than diethyl ether; however, you may need to determine the appropriate solvent by experimentation. Once you have selected a solvent, this same solvent should be used throughout the procedure when an organic solvent is required. If you use diethyl ether, you must use two steps to dry the organic layer. First, the organic layer must be mixed with saturated sodium chloride (see Technique 12, Section 12.9, Saturated Salt Solution), and then the liquid is dried over granular anhydrous sodium sulfate. (see Technique 12, Section 12.9, Drying Procedure with Anhydrous Sodium Sulfate).

¹ Your mixture may be one of the following: (1) 50% benzoic acid, 40% benzoin, 10% 1,4-dibromobenzene; (2) 50% fluorene, 40% *o*-toluic acid, 10% 1,4-dibromobenzene; (3) 50% phenanthrene, 40% methyl 4-aminobenzoate, 10% 1,4-dibromobenzene; or (4) 50% 4-aminoacetophenone, 40% 1,2,4,5-tetrachlorobenzene, 10% 1,4-dibromobenzene. Other mixtures are given in the Instructor's Manual, along with some suggestions about these mixtures.

Purification

To improve the purity of your final samples, it may be helpful to include a backwashing step at the appropriate place in your procedure. See Technique 12, Section 12.11, for a discussion of this technique. Crystallization will likely be required to purify both of the compounds you isolate. To find an appropriate solvent, you should consult a handbook. You can also use the procedure in Technique 12, Section 11.6, to determine a good solvent experimentally. Note that diethyl ether or other very low boiling solvents are not generally good solvents for performing crystallizations. If you use water as a solvent, you will need to let the crystals air-dry overnight. Take melting points of your final samples to determine if you have obtained both compounds in a pure form. Hand in each compound in a labeled vial.

When performing the laboratory work, you should strive to obtain a high recovery of both compounds in a highly pure form. If your procedure fails, modify it and repeat the experiment.

REPORT

Write out a complete procedure by which you separated and isolated pure samples of two of the compounds in your mixture. Describe how you determined that your procedure was successful and give any data or results used for this purpose. Calculate the percentage recovery for both compounds.

You should include your outline of the separation scheme with your report (see Technique 12, Section 12.12).

Chromatography

Thin-layer chromatography

Column chromatography

Following a reaction with thin-layer chromatography

Chromatography is perhaps the most important technique used by organic chemists to separate the components of a mixture. This technique involves the distribution of the different compounds or ions in the mixture between two phases, one of which is stationary and the other moving. Chromatography works on much the same principle as solvent extraction. In extraction, the components of a mixture are distributed between two solvents according to their relative solubilities in the two solvents. The separation process in chromatography depends on differences in how strongly the components of the mixture are adsorbed to the stationary phase and how soluble they are in the moving phase. These differences depend primarily on the relative polarities of the components in the mixture.

There are many types of chromatographic techniques, ranging from thin-layer chromatography, which is relatively simple and inexpensive, to high-performance liquid chromatography, which is very sophisticated and expensive. In this experiment, you will use two of the most widely used chromatographic techniques: thin-layer and column chromatography. The purpose of this experiment is to give you practice in performing these two techniques, to illustrate the principles of chromatographic separations, and to demonstrate how thin-layer and column chromatography are used in organic chemistry.

REQUIRED READING

New: Technique 19 Column Chromatography
Technique 20 Thin-Layer Chromatography

SPECIAL INSTRUCTIONS

Many flammable solvents are used in this experiment. Use Bunsen burners for making micropipettes in a part of the lab that is separate from where the solvents are being used. The thin-layer chromatography should be performed in the hood.

SUGGESTED WASTE DISPOSAL

Dispose of methylene chloride in the container designated for halogenated organic wastes. Dispose of all other organic solvents in the container for nonhalogenated organic solvents. Place the alumina in the container designated for wet alumina.

NOTES TO THE INSTRUCTOR

The column chromatography should be performed with activated alumina from EM Science (No. AX0612-1). The particle sizes are 80–200 mesh, and the material is Type F-20. The alumina should be dried overnight in an oven at 110°C and stored in a tightly sealed bottle. Alumina more than several years old may need to be dried for a longer time at a higher temperature.


For thin-layer chromatography (TLC), use flexible silica gel plates from Whatman with a fluorescent indicator (No. 4410 222). If the TLC plates have not been purchased recently, they should be placed in an oven at 100°C for 30 minutes and stored in a desiccator until used. If you use different alumina or different thin-layer plates, try out the experiment before using it with a class. Other materials than those specified here may give different results from those indicated in this experiment.

Grind up the fluorenone flakes into smaller pieces for easier dispensing. Commercially available fluorenone is often contaminated with fluorene and fluorene, and fluorenone is often contaminated with fluorene. If iodine is used to visualize the spots in Experiment 6A, these contaminants will likely be invisible. However, if a UV lamp, which is more sensitive, is used, the contaminants will likely be visible. These compounds can be purified by crystallization (see Instructor's Manual) and then the contaminants will likely be invisible even when the spots are visualized under a UV lamp. It is best to use iodine to visualize the spots in Experiment 6C even if the fluorenone is pure. Since iodine is not as sensitive as a UV lamp, students will observe a more gradual change in the intensities of the spots for the two compounds when iodine is used.

6A EXPERIMENT 6A

Thin-Layer Chromatography

In this experiment, you will use thin-layer chromatography (TLC) to separate a mixture of three compounds: fluorene, fluorenone, and fluorenone:

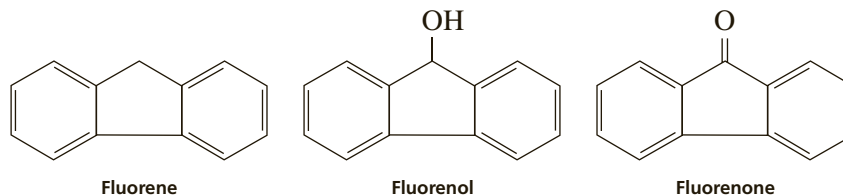


Based on the results with known samples of these compounds, you will determine which compounds are found in an unknown sample. Using TLC to identify the components in a sample is a common application of this technique.

PROCEDURE

Preparing the TLC Plate

Technique 20 describes the procedures used for thin-layer chromatography. Use a 10 cm × 5.3 cm TLC plate (Whatman Silica Gel Plates No. 4410 222). These plates have a flexible backing but should not be bent excessively. They should be handled carefully

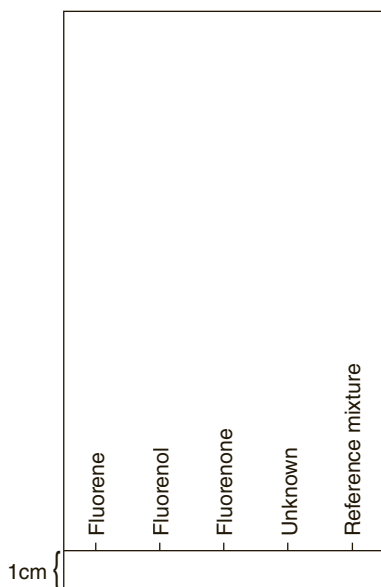

Preparation of a TLC plate.

or the adsorbent may flake off them. Also, they should be handled only by the edges; the surface should not be touched. Using a lead pencil (not a pen), lightly draw a line across the plate (short dimension) about 1 cm from the bottom (see figure). Using a centimeter ruler, move its index about 0.6 cm in from the edge of the plate and lightly mark off five 1-cm intervals on the line. These are the points at which the samples will be spotted.

Prepare five micropipettes to spot the plate. The preparation of these pipettes is described and illustrated in Technique 20, Section 20.4. Prepare a TLC development chamber with methylene chloride (see Technique 20, Section 20.5). A beaker covered with aluminum foil or a wide-mouth, screw-cap bottle is a suitable container to use (see Technique 20, Figure 20.4). The backing on the TLC plates is thin, so if it touches the filter paper liner of the development chamber at any point, solvent will begin to diffuse onto the adsorbent surface at that point. To avoid this, be sure that the filter paper liner does not go completely around the inside of the container. A space about 2.5 inches wide must be provided. (*Note:* This development chamber will also be used for Experiments 6C and 6D.)

On the plate, starting from left to right, spot fluorene, fluorenol, fluorenone, the unknown mixture, and the standard reference mixture, which contains all three compounds.¹ For each of the five samples, use a

different micropipette to spot the sample on the plate. The correct method of spotting a TLC plate is described in Technique 20, Section 20.4. Take up part of the sample in the pipette (don't use a bulb; capillary action will draw up the liquid). Apply the sample by touching the pipette *lightly* to the thin-layer plate. The spot should be no larger than 2 mm in diameter. It will usually be sufficient to spot each sample once or twice.² If you need to spot the sample more than once, allow the solvent to evaporate completely between successive applications and spot the plate in exactly the same position each time. Save the samples in case you need to repeat the TLC.³

Developing the TLC Plate

Place the TLC plate in the development chamber, making sure that the plate does not come in contact with the filter paper liner. Remove the plate when the solvent front is 1–2 cm from the top of the plate. Using a lead pencil, mark the position of the solvent front. Set the plate on a piece of paper towel to dry. When the plate is dry, place the plate in a jar containing a few iodine crystals, (see Technique 20, Section 20.7) cap the jar, and leave it in the jar until the spots begin to appear. Remove the plate from the jar and lightly outline all the spots that became visible with the iodine treatment. Using a ruler marked in millimeters, measure the distance that each spot has traveled relative to the solvent front. Calculate the R_f values for each spot (see Technique 20, Section 20.9). Explain the relative positions of the three compounds in terms of their polarities. Identify the compound or compounds that are found in the unknown mixture. At the instructor's option, submit the TLC plate with your report.

¹ *Note to the instructor:* The individual compounds and the reference mixture containing all three compounds are prepared as 2% solutions in acetone. The unknown mixture may contain one, two, or all three of the compounds dissolved in acetone.

² If a UV lamp will be used to visualize the spots after developing the plate, you should spot each sample only once.

³ After you have developed the plate and seen the spots, you will be able to tell if you need to rerun the TLC plate. If the spots are too faint to see clearly, you need to spot the sample more. If any of the spots show tailing (Technique 19, Section 19.12), then less sample is needed.

6B EXPERIMENT 6B

Selecting the Correct Solvent for Thin-Layer Chromatography

In Experiment 6A, you were told what solvent to use for developing the TLC plate. In some experiments, however, it will be necessary to determine an appropriate development solvent by experimentation (Technique 20, Section 20.6). In this experiment, you will be instructed to try three solvents for separating a pair of related compounds that differ slightly in polarity. Only one of these solvents will separate the two compounds enough so that they can be easily identified. For the other two solvents, you will be asked to explain, in terms of their polarities, why they failed.

PROCEDURE

Preparation

Your instructor will assign you a pair of compounds to run on TLC, or you will select your own pair.⁴ You will need to obtain about 0.5 mL of three solutions: one solution of each of the two individual compounds and a solution containing both compounds. Prepare three thin-layer plates in the same way as you did in Experiment 5A, except that each plate should be 10 cm × 3.3 cm. When you mark them with a pencil for spotting, make three marks 1 cm apart. Prepare three micropipettes to spot the plates. Prepare three TLC development chambers as you did in Experiment 5A, with each chamber containing one of the three solvents suggested for your pair of compounds.

Developing the TLC Plate

On each plate, spot the two individual compounds and the mixture of both compounds. For each of the three samples, use a different micropipette to spot the sample on the plates. Place each TLC plate in one of the three development chambers, making sure that the plate does not come in contact with the filter paper liner. Remove each plate when the solvent front is 1–2 cm from the top of the plate. Using a lead pencil, mark the position of the solvent front. Set the plate on a piece of paper towel to dry. When the plate is dry, observe it under a short-wavelength UV lamp, preferably in a darkened hood or a darkened room. With a pencil, lightly outline any spots that appear. Next, place the plate in a jar containing a few iodine crystals, and leave it in the jar until the spots begin to appear. Remove the plate from the jar and lightly outline all the spots that became visible with the iodine treatment. Using a ruler marked in millimeters, measure the distance that each spot has traveled relative to the solvent front. Calculate the R_f values for each spot. At the instructor's option, submit the TLC plates with your report.


Which of the three solvents resolved the two compounds successfully? For the two solvents that did not work, explain, in terms of their polarities, why they failed.

⁴ *Note to the instructor:* Possible pairs of compounds are given in the following list. The two compounds to be resolved are given first, followed by the three developing solvents to try: (1) benzoin and benzil; acetone, methylene chloride, hexane; (2) vanillin and vanillyl alcohol; acetone, 50% toluene–50% ethyl acetate, hexane; (3) diphenylmethanol and benzophenone; acetone, 70% hexane–30% acetone, hexane. Each compound in a pair should be prepared individually and as a mixture of the two compounds. Prepare all of them as 1% solutions in acetone.

6C EXPERIMENT 6C

Monitoring a Reaction with Thin-Layer Chromatography

Thin-layer chromatography is a convenient method for monitoring the progress of a reaction (Technique 20, Section 20.10). This technique is especially useful when the appropriate reaction conditions have not yet been worked out. By using TLC to follow the disappearance of a reactant and the appearance of a product, it is relatively easy to decide when the reaction is complete. In this experiment, you will monitor the reduction of fluorenone to fluorenlol:

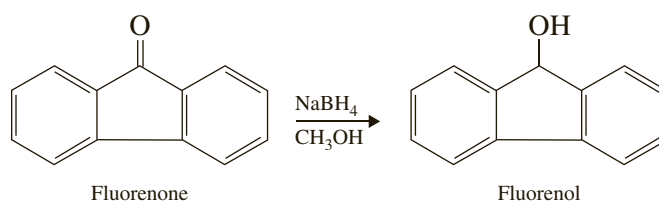

Although the appropriate reaction conditions for this reaction are already known, using TLC to monitor the reaction will demonstrate how to use this technique.

PROCEDURE

Preparation

Work with a partner on this part of the experiment. Prepare two thin-layer plates in the same way as you did in Experiment 6A, except that one plate should be 10 cm × 5.3 cm and the other one, 10 cm × 4.3 cm. When you mark them with a pencil for spotting, make five marks 1 cm apart on the first plate and four marks on the second plate. During the reaction, you will be taking five samples from the reaction mixture at 0, 15, 30, 60, and 120 seconds. Three of these samples should be spotted on the larger plate and two of them on the smaller one. In addition, each plate should be spotted with two reference solutions, one containing fluorenone and the other fluorenlol. Using a pencil to make very light marks, indicate at the top of each plate where each sample will be spotted so that you can keep track of them. Write the number of seconds and an abbreviation for the two reference compounds. Use the same TLC development chamber with methylene chloride that you used in Experiment 5A. Prepare seven micropipettes to spot the plates.

Running the Reaction

Once sodium borohydride has been added to the reaction mixture (see next paragraph), take samples at the times just indicated. Because this must be done in such a short time, you must be well prepared before starting the reaction. One person should be the timekeeper and the other person should take the samples and spot the plates. Spot each sample once, using a different pipette for each sample.

Place a magnetic spin vane (Technique 7, Figure 7.4A) into a 5-mL conical vial. Add 0.20 g of fluorenone and 4 mL of methanol to the vial. Place the vial on a magnetic stirrer, using either an aluminum block or a clamp to hold the vial in

place. Stir the mixture until all the solid has dissolved. Now take the first sample (the “0 second” sample) and spot the plate. Using smooth weighing paper, weigh 0.020 g of sodium borohydride⁵ and immediately add it to the reaction mixture. If you wait too long to add it, the sodium borohydride will become sticky because it absorbs moisture from the air. Begin timing the reaction as soon as the sodium borohydride is added. Use the micropipettes to remove samples of the reaction mixture at the following times: 15, 30, 60, and 120 seconds. Use a different micropipette each time and spot a TLC plate with each sample. On each plate, also spot the two reference solutions of fluorenone and fluorenol in acetone. After developing the plates and allowing them to dry, visualize the spots with iodine, as described in Part A. Make a sketch of your plates and record the results in your notebook. Do these results indicate that the reaction went to completion? In addition to the TLC results, what other visible evidence indicated that the reaction went to completion? Explain.

Isolation of Fluorenol (optional procedure)

Using a Pasteur pipette, transfer the reaction mixture to a 10-mL Erlenmeyer flask. Add 1 mL of water and heat the mixture almost to boiling for about 2 minutes. Allow the flask to cool slowly to room temperature in order to crystallize the product. Then place the flask in an ice-water bath for several minutes to complete crystallization. Collect the crystals by vacuum filtration, using a small Hirsch funnel (Technique 8, Section 8.3). Wash the crystals with three 1.0-mL portions of an ice-cold mixture of 80% methanol and 20% water. After the crystals are dry, weigh them and determine their melting point (literature, 153–154°C).

6D EXPERIMENT 6D

Column Chromatography

The principles of column chromatography are similar to those of thin-layer chromatography. The primary difference is that the moving phase in column chromatography travels downward, whereas in TLC the solvent ascends the plate. Column chromatography is used more often than TLC to separate relatively large amounts of compounds. With column chromatography, it is possible to collect pure samples of the separated compounds and perform additional tests on them.

In this experiment, fluorene and fluorenone will be separated by column chromatography using alumina as the adsorbent. Because fluorenone is more polar than fluorene, fluorenone will be absorbed to the alumina more strongly. Fluorene will elute off the column with a nonpolar solvent hexane, whereas fluorenone will not come off until a more polar solvent (30% acetone–70% hexane) is put on the column. The purities of the two separated compounds will be tested by TLC and melting points.


⁵ *Note to the Instructor:* The sodium borohydride should be checked to see whether it is active: Place a small amount of powdered material in some methanol and heat it gently. If the hydride is active, the solution should bubble vigorously. If using an old bottle, it is also good to check the material for stickiness due to absorption of water. If it is too sticky, it can be difficult for students to weigh it out.

PROCEDURE

Advance Preparation

Before running the column, assemble the following glassware and liquids. Obtain four dry test tubes (16 × 100 mm) and number them 1 through 4. Prepare two dry Pasteur pipettes with bulbs attached. Place 9.0 mL of hexane, 2.0 mL of acetone, and 2.0 mL of a solution of 70% hexane–30% acetone (by volume) into three Erlenmeyer flasks. Clearly label and stopper each flask. Place 0.3 mL of a solution containing fluorene and fluorenone into a small test tube.⁶ Stopper the test tube. Prepare one 10 cm × 3.3 cm TLC plate with four marks for spotting. Use the same TLC development chamber with methylene chloride that you used in Part A. Prepare four micropipettes to spot the plates.

Prepare a chromatography column packed with alumina. Place a loose plug of cotton in a Pasteur pipette (5¾-inch) and push it gently into position using a glass rod (see figure for the correct position of the cotton). *Do not ram the cotton tightly, because this may result in the solvent flowing through the column too slowly.* Using a file, score the Pasteur pipette about 1 cm below the cotton plug. To break the tip off the pipette, put your thumbs together at the place on the pipette that you scored and push quickly with both thumbs.

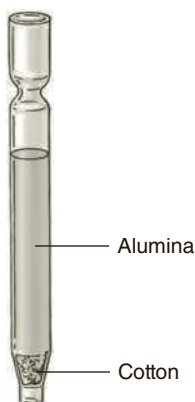

Chromatography column.

CAUTION



Wear gloves or use a towel to protect your hands from being cut while breaking the pipette.

Add 1.25 g of alumina (EM Science, No. AX0612-1) to the pipette while tapping the column gently with your finger.⁷ When all the alumina has been added, tap the column with your finger for several seconds to ensure that the alumina is tightly packed. Clamp the column in a vertical position so that the bottom of the column is just above the height of the test tubes you will be using to collect the fractions. Place test tube 1 under the column.

Running the Column

Using a Pasteur pipette, add 3 mL of hexane to the column. The column must be completely moistened by the solvent. Drain the excess hexane until the level of hexane reaches the top of the alumina. Once hexane has been added to the alumina, the top of the column must not be allowed to run dry. If necessary, add more hexane.

NOTE: It is essential that the liquid level not be allowed to drain below the surface of the alumina at any point in this procedure.

⁶ *Note to the instructor:* This solution should be prepared for the entire class by dissolving 0.3 g of fluorene and 0.3 g of fluorenone in 9.0 mL of a mixture of 5% methylene chloride–95% hexane. Store this solution in a closed container to prevent evaporation of solvent. This will provide enough solution for 20 students, assuming little spillage or other types of waste.

⁷ As an option, students may prepare a microfunnel from a 1-mL disposable plastic pipette. The microfunnel is prepared by (1) cutting the bulb in half with a scissors and (2) cutting the stem at an angle about ½ inch below the bulb. This funnel can be placed in the top of the column (Pasteur pipette) to aid in filling the column with alumina or with the solvents (see Technique 19, Section 19.6).

When the level of the hexane reaches the top of the alumina, add the solution of fluorene and fluorenone to the column using a Pasteur pipette. Begin collecting the eluent in test tube 2. Just as the solution penetrates the column, add 1 mL of hexane and drain until the surface of the liquid has reached the alumina. Add another 5 mL of hexane. As fluorene elutes off the column, some solvent will evaporate, leaving solid fluorene on the tip of the pipette. Using a Pasteur pipette, dissolve this solid off the column with a few drops of acetone. It may be necessary to do this several times, and the acetone solution is also collected in tube 2.

After you have added all the hexane, change to the more polar solvent (70% hexane–30% acetone).⁸ When changing solvents, do not add the new solvent until the last solvent has nearly penetrated the alumina. The yellow band (fluorenone) should now move down the column. Just before the yellow band reaches the bottom of the column, place test tube 3 under the column. When the eluent becomes colorless again, place test tube 4 under the column and stop the procedure.

Tube 2 should contain fluorene and tube 3, fluorenone. Test the purities of these two samples using TLC. You must spot the solution from tube 2 several times in order to apply enough sample on the plate to be able to see the spots. On the plate, also spot a reference solution containing fluorene and fluorenone. After developing the plate and allowing it to dry, visualize the spots with iodine. What do the TLC results indicate about the purities of the two samples?

Using a warm water bath (40–60°C) and a stream of nitrogen gas or air, evaporate the solvent from test tubes 2 and 3. As soon as all the solvent has evaporated from each of the tubes, remove them from the water bath. There may be a yellow oil in tube 3, but it should solidify when the tube cools to room temperature. If it does not, cool the tube in an ice-water bath and scratch the bottom of the test tube with a glass stirring rod or a spatula. Determine the melting points of the fluorene and fluorenone. The melting point of fluorene is 116–117°C and of fluorenone is 82–85°C.

REPORT

Experiment 6A

1. Calculate the R_f values for each spot. Include the actual plate or a sketch of the plate with your report.
2. Explain the relative R_f values for fluorene, fluorenone, and fluorenone in terms of their polarities and structures.
3. Give the composition of the unknown that you were assigned.

Experiment 6B

1. Record the names and structures of the two compounds that you ran on TLC.
2. Which solvent resolved the two compounds successfully?
3. For the other two solvents, explain, in terms of their polarities, why they failed.

Experiment 6C

1. Make a sketch of the TLC plate or include the actual plate with your report. Interpret the results. When was the reaction complete?
2. What other visible evidence indicated that the reaction went to completion?
3. If you isolated the fluorenone, record the melting point and the weight of this product.

⁸ Sometimes the fluorenone also moves through the column with hexane. Therefore, be sure to change to test tube 3 if the yellow band starts to emerge from the column.

Experiment 6D

1. Describe the TLC results on the samples in test tubes 2 and 3. What does this indicate about the purities of the two samples?
2. Record the melting points for the dried solids found in tubes 2 and 3. What do they indicate about the purities of the two samples?

QUESTIONS

1. Each of the solvents given should effectively separate one of the following mixtures by TLC. Match the appropriate solvent with the mixture that you would expect to separate well with that solvent. Select your solvent from the following: hexane, methylene chloride, or acetone. You may need to look up the structures of solvents and compounds in a handbook.
 - a. 2-Phenylethanol and acetophenone
 - b. Bromobenzene and *p*-xylene
 - c. Benzoic acid, 2,4-dinitrobenzoic acid, and 2,4,6-trinitrobenzoic acid
2. The following questions relate to the column chromatography experiment performed in Experiment 6D.
 - a. Why does the fluorene elute first from the column?
 - b. Why was the solvent changed in the middle of the column procedure?
3. Consider the following errors that could be made when running TLC. Indicate what should be done to correct the error.
 - a. A two-component mixture containing 1-octene and 1,4-dimethylbenzene gave only one spot with an R_f value of 0.95. The solvent used was acetone.
 - b. A two-component mixture containing a dicarboxylic acid and tricarboxylic acid gave only one spot with an R_f value of 0.05. The solvent used was hexane.
 - c. When a TLC plate was developed, the solvent front ran off the top of the plate.

Infrared Spectroscopy and Boiling-Point Determination

Infrared spectroscopy

Boiling-point determination

Organic nomenclature

Critical thinking application

The ability to identify organic compounds is an important skill that is frequently used in the organic laboratory. Although there are several spectroscopic methods and many chemical and physical tests that can be used for identification, the goal of this experiment is to identify an unknown liquid using infrared spectroscopy and a boiling-point determination. Both methods are introduced in this experiment.

REQUIRED READING

<i>New:</i> Technique 4	How to Find Data for Compounds: Handbooks and Catalogs
Technique 13	Physical Constants of Liquids: The Boiling Point and Density, Part A. "Boiling Points and Thermometer Correction"
Technique 25	Infrared Spectroscopy

SPECIAL INSTRUCTIONS

Many of the unknown liquids used for this experiment are flammable; therefore, do not use any flames in the laboratory. Also, be careful when handling all of the liquids because many of them are potentially toxic.

This experiment can be performed individually, with each student working on one unknown. However, the opportunity to learn is greater if students work in groups of three. In this case, each group is assigned three unknowns. Each student in the group obtains an infrared spectrum and performs a boiling-point determination on one of the unknowns. Subsequently, the student shares this information with the other two students in the group. Then each student analyzes the collective results for the three unknowns and writes a laboratory report based on all three unknowns. Your instructor will inform you whether you should work alone or in groups.

SUGGESTED WASTE DISPOSAL

If you have not identified the unknown by the end of the laboratory period, you should return the unknown liquid to your instructor in the original container in which it was issued to you. If you have identified the compound, dispose of it in either the container for halogenated waste or the one for nonhalogenated waste, whichever is appropriate.

NOTES TO THE INSTRUCTOR

If you choose to have students work in groups of three, be sure to assign unknowns that differ both in structure *and* functional group, with at least one aromatic compound in each set. If the experiment is performed early in the year, students may have some difficulty in drawing the structures of the compounds that are in the list of possible unknowns, and they will need help. For each unknown, structures will be needed for several of the possible compounds. In fact, compounds with boiling points as much as 5°C higher than the experimental boiling point should be considered because student-determined boiling points are frequently low. This will depend on the method used and the skill of the person performing the technique. *The Merck Index*, the *CRC Handbook of Chemistry and Physics*, and the lecture textbook can all be helpful in determining these structures. Technique 4, "How to Find Data for Compounds: Handbooks and Catalogs," provides helpful information for students just beginning to use handbooks. The nuclear magnetic resonance (NMR) portion of the experiment is optional. We suggest that access to the NMR be granted only after a plausible solution has been tendered. If you do not have an NMR, there are several on-line databases where you can obtain a printed copy of the spectrum to hand to students.

For the boiling point determination, we prefer the semimicroscale direct method described in Technique 13, Section 13.2, Semimicroscale Direct Method. A convenient, safe, and accurate way to determine the boiling point is to use a Vernier LabQuest with a stainless steel probe or a Vernier LabPro interface with a laptop computer and stainless steel probe (see Technique 13, Section 13.4, and Figure 13.7). With both of these methods, students observe a graph of time vs. temperature. Another convenient method is to use a digital thermometer with a stainless steel probe; however, most of these devices do not provide a graph of the temperature (see Technique 13, Section 13.4, and Figure 13.7). All of these methods are more accurate than using non-mercury thermometers. If a glass thermometer is used, the boiling point will be most accurate if a partial immersion mercury thermometer is used. See the Instructor's Manual for additional comments about the use of these different methods.

PROCEDURE

Part A. Infrared Spectrum

Obtain the infrared spectrum of your unknown liquid (Technique 25, Section 25.2). If you are working in a group, provide copies of your spectrum for everyone in your group. Identify the significant absorption peaks by labeling them *right on the*

spectrum, and include the spectrum in your laboratory report. Absorption peaks corresponding to the following groups should be identified:

- C—H (sp^3)
- C—H (sp^2)
- C—H (aldehyde)
- O—H
- C=O
- C=C (aromatic)
- aromatic substitution pattern
- C—O
- C—X (if applicable)
- N—H

Part B. Boiling-Point Determination

Perform a boiling-point determination on your unknown liquid (Technique 13, Section 13.2). Your instructor will indicate which method to use. Depending on the method used and the skill of the person performing the technique, boiling points can sometimes be slightly inaccurate. When experimental boiling points are inaccurate, it is most common for them to be lower than the literature value. The difference may be as much as 5°C , especially for higher-boiling liquids and if you use a non-mercury thermometer. If you use a digital thermometer with a stainless steel probe or one of the Vernier devices with a stainless steel probe, the results should be within $2\text{--}3^\circ\text{C}$ of the actual boiling point. Results obtained using a partial immersion mercury thermometer also tend to be very good. Your instructor may be able to give you more guidance about what level of accuracy you can expect.

Part C. Analysis and Report

Using the structural information from the infrared spectrum and the boiling point of your unknown, identify this liquid from the list of compounds in the table. If you are working in a group, you will need to do this for all three compounds. In order to make use of the structural information determined from the infrared spectrum, you will need to know the structures of the compounds that have boiling points close to the value you experimentally determined. You may need to consult *The Merck Index* or the *CRC Handbook of Chemistry and Physics*. It may also be helpful to look up these compounds in the index of your lecture textbook. If there is more than one compound that fits the infrared spectrum and is within a few degrees of the experimental boiling point, you should list all of these in your laboratory report.

In your laboratory report, include (1) the infrared spectrum with the significant absorption peaks identified *right on the spectrum*, (2) the experimental boiling point for your unknown, and (3) your identification of the unknown. Explain your justifications for making this identification and write out the structure of this compound.

Optional Exercise: NMR Spectrum

At the option of your instructor, you may be asked to determine the nuclear magnetic resonance spectrum of your unknown liquid (Technique 26, Section 26.1). Alternatively, your instructor may issue you a previously run spectrum of your compound. You should provide structural assignments for all of the groups of hydrogens that are present. Do this *right on the spectrum*. If you have correctly determined the identity of your unknown, all the groups of hydrogens (and their chemical shifts) should fit your structure. Include the properly labeled spectrum in your report and explain why it fits the suggested structure.

List of possible unknown liquids

Compound	BP (°C)	Compound	BP (°C)
Acetone	56	Butyl acetate	127
2-Methylpentane	62	2-Hexanone	128
sec-Butylamine	63	Morpholine	129
Isobutyraldehyde	64	3-Methyl-1-butanol	130
Methanol	65	Hexanal	130
Isobutylamine	69	Chlorobenzene	132
Hexane	69	2,4-Pentanedione	134
Vinyl acetate	72	Cyclohexylamine	135
1,3,5-Trifluorobenzene	75	Ethylbenzene	136
Butanal	75	p-Xylene	138
Ethyl acetate	77	1-Pentanol	138
Butylamine	78	Propionic acid	141
Ethanol	78	Pentyl acetate	142
2-Butanone	80	4-Heptanone	144
Cyclohexane	81	2-Ethyl-1-butanol	146
Isopropyl alcohol	82	N-Methylcyclohexylamine	148
Cyclohexene	83	2,2,2-Trichloroethanol	151
Isopropyl acetate	85	2-Heptanone	151
Triethylamine	89	Heptanal	153
3-Methylbutanal	92	Isobutyric acid	154
3-Methyl-2-butanone	94	Bromobenzene	156
1-Propanol	97	Cyclohexanone	156
Heptane	98	Dibutylamine	159
tert-Butyl acetate	98	Cyclohexanol	160
2,2,4-Trimethylpentane	99	Butyric acid	162
2-Butanol	99	Furfural	162
Formic acid	101	Diisobutyl ketone	168
2-Pentanone	101	Furfuryl alcohol	170
2-Methyl-2-butanol	102	Octanal	171
Pentanal	102	Decane	174
3-Pentanone	102	Isovaleric acid	176
Propyl acetate	102	Limonene	176
Piperidine	106	1-Heptanol	176
2-Methyl-1-propanol	108	Benzaldehyde	179
1-Methylcyclohexene	110	Cycloheptanone	181
Toluene	111	1,4-Diethylbenzene	184
sec-Butyl acetate	111	Iodobenzene	186
Pyridine	115	1-Octanol	195
4-Methyl-2-pentanone	117	Methyl benzoate	199
2-Ethylbutanal	117	Methyl phenyl ketone	202
Methyl 3-methylbutanoate	117	Benzyl alcohol	204
Acetic acid	118	4-Methylbenzaldehyde	204
1-Butanol	118	Ethyl benzoate	212
Octane	126		

Simple and Fractional Distillation

Simple distillation

Fractional distillation

Gas chromatography

Distillation is a technique frequently used to separate and purify a liquid component from a mixture. Simply stated, distillation involves heating a liquid mixture to its boiling point, where liquid is rapidly converted to vapor. The vapors, richer in the more volatile component, are then condensed into a separate container. When the components in the mixture have sufficiently different vapor pressures (or boiling points), they can be separated by distillation.

The purpose of this experiment is to illustrate the use of distillation for separating a mixture of two volatile liquids with different boiling points. Each mixture, which will be issued as an unknown, will consist of two liquids from the following table.

Compound	Boiling Point (°C)
Hexane	69
Cyclohexane	80.7
Heptane	98.4
Toluene	110.6

The liquids in the mixture will be separated by two distillation techniques: simple and fractional distillation. The results of these two methods will be compared by analyzing the composition of the **distillate** (the distilled liquid) using gas chromatography. You will also construct a graph of the distillation temperature versus the total volume of distillate collected. This graph will allow you to determine the approximate boiling points of the two liquids and to make a graphical comparison of the two distillation methods.

Experiment 8A is designed to be performed with semimicroscale glassware using a conventional distillation apparatus. A microscale alternative with a Hickman head is given in Experiment 8B; however, the scale is the same in both cases, and the experiment can be performed more easily with the semimicroscale glassware.

REQUIRED READING

<i>New:</i> Technique 14	Simple Distillation
Technique 15	Fractional Distillation
Technique 22	Gas Chromatography

¹ This experiment is based on a similar one developed by James Patterson, North Seattle Community College, Seattle.

SPECIAL INSTRUCTIONS

Many flammable solvents are used in this experiment; therefore, do not use any flames in the laboratory.

Work in pairs on this experiment. Each pair of students will be assigned an unknown containing two liquids found in the table above. One student in the pair should perform a simple distillation and the other student, a fractional distillation. The results from these two methods will be compared.

SUGGESTED WASTE DISPOSAL

Dispose of all organic liquids in the container for nonhalogenated organic solvents.

NOTES TO THE INSTRUCTOR

One method of insulating the air condenser used for the fractional distillation column is provided by employing two layers of clear flexible tubing (PVC, polyvinyl chloride) over the air condenser. For a ½-in. diameter column, use ½-in. × ⅝-in. outer-diameter plastic tubing on the inside and ⅝-in. I.D. × ⅞-in. O.D. tubing on the outside. Cut the tubing into 3½-in. lengths. Make a slit from end to end so that the lengths can slip over the column. Slit the tubing using sharp scissors or a razor knife with a proper handle. Do not use a razor blade, or you may get badly cut. The clear tubing lets you see what is going on in the column and also provides some insulation. Another method of insulating the fractionating column is to wrap the air condenser with a cotton pad about 3½ in. square. Prepare the cotton pad by covering both sides of one layer of cotton with aluminum foil. Wrap this pad entirely with duct tape to hold the cotton in place and to make a more durable pad. Hold the pad in place with tape or twist ties.

A convenient, safe, and accurate way to monitor the temperature during the distillation is to use a Vernier LabQuest device with a stainless steel probe or a Vernier LabPro interface with a laptop computer and stainless steel probe (see Technique 13, Section 13.4 and Technique 14, Section 14.5 and Figure 14.12). With both of these methods, students observe a graph of time vs. temperature. Another convenient method is to use a digital thermometer with a stainless steel probe (see Technique 14, Section 14.5 and Figure 14.12); however, most of these devices do not provide a graph of the temperature. All of these methods are more accurate than using non-mercury thermometers. If a glass thermometer is used, the temperature will be most accurate if a partial immersion mercury thermometer is used. See the Instructor's Manual for additional comments about the use of these devices, including suitable stainless steel probes for this experiment.

Prepare unknown mixtures consisting of the following pairs of liquids: hexane–heptane, hexane–toluene, and cyclohexane–toluene. Distillation of these mixtures should provide a good contrast between the two distillation methods. See the Instructor's Manual for directions on how to make up these mixtures.

Unless the samples are analyzed by gas chromatography immediately after the distillation, it is essential that the samples be stored in leak-proof vials. We have found that GC-MS vials work much better for this purpose than conical vials found in microscale glassware kits.

The gas chromatograph is prepared as follows: column temperature, 140°C; injection temperature, 150°C; detector temperature, 140°C; carrier gas flow rate,

100 mL/min. The recommended column is 8 feet long, with a stationary phase such as Carbowax 20M.

You should determine response factors for the four liquids given in the table provided at the beginning of this experiment. Because the data in this experiment are expressed as volume, the response factors should also be based on volume. Inject a mixture containing equal volumes of all four compounds and determine the relative peak areas. Choose one compound as the standard and define its response factor to be equal to 1.00. Calculate the other response factors based on this reference. Typical response factors are given in Footnote 4. Although this experiment can be performed with a Hickman head (see Experiment 8B), it is more difficult to control the distillation rate, and the results may not be as satisfactory as those obtained in Experiment 8A.

8A EXPERIMENT 8A

Simple and Fractional Distillation (Semimicroscale Procedure)

PROCEDURE

You should work in pairs on this experiment. Each pair of students will be assigned an unknown mixture containing equal volumes of two of the liquids from the table at the beginning of this experiment. One student should perform a simple distillation on the mixture, and the other student should perform a fractional distillation.

Apparatus


The temperature during the distillation may be monitored either with a thermometer or a stainless steel temperature probe. If a stainless steel probe is used, it must be used in conjunction with either a digital thermometer or one of the Vernier devices (see Technique 13, Section 13.4, and Technique 14, Figure 14.12). Your instructor will provide instructions about the method that you will use. Assemble the appropriate distillation apparatus (see figures). Carefully notice the position of the thermometer in these figures. The bulb of the thermometer or the end of the stainless steel probe must be placed well below the sidearm or they will not read the temperature correctly. If a thermometer is used, it is held in place with a thermometer adapter. If a temperature probe is used, it is held in place with a rubber septum (see Technique 14, Figure 14.12).

If performing the fractional distillation, pack the air condenser uniformly with 0.8–0.9 g of stainless steel cleaning pad material. The easiest way to pack the column is to cut several strands of the cleaning pad with the correct weight. Using a long wire with a bend end, pull the cleaning pad through the condenser. After releasing the wire, use a metal spatula or glass stirring rod to adjust the position of the cleaning pad. Do not pack the material too tightly at any one place in the condenser.

CAUTION



You should wear heavy cotton gloves when handling the stainless steel cleaning pad. The edges are very sharp and can easily cut into the skin.

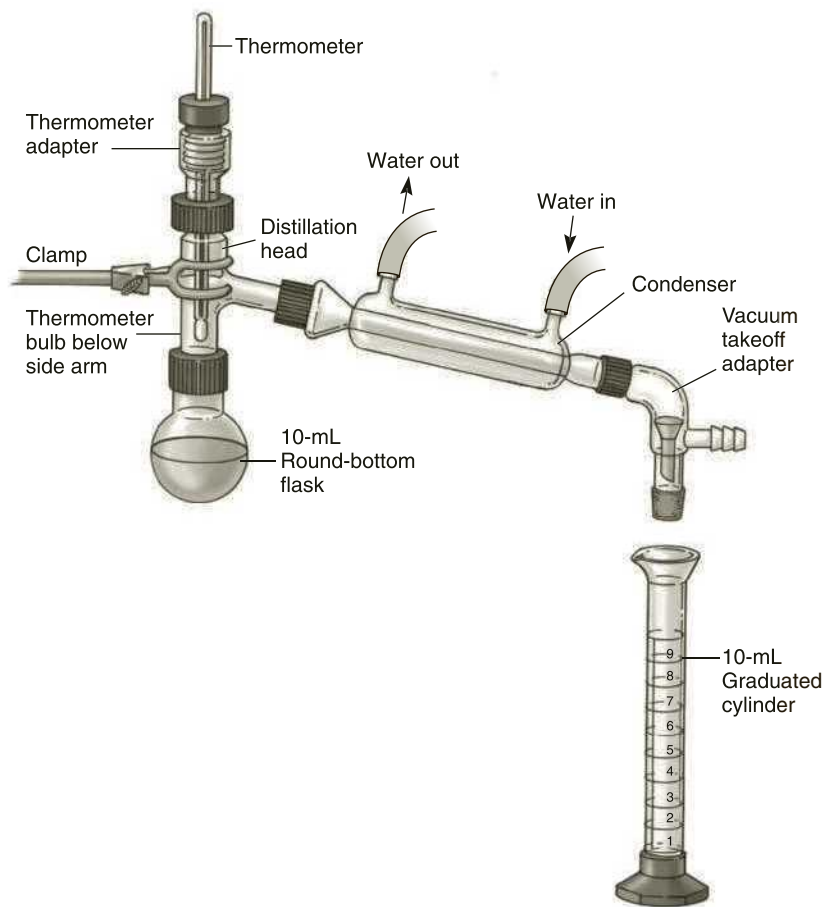


Apparatus for simple distillation. A digital thermometer can also be used; see Section 13.4 and Figure 14.12 in the Techniques.

Wrap the glass section of the air condenser between the two plastic caps with plastic tubing as described in Notes to the Instructor. Alternatively, use the method with a cotton pad (see Notes to the Instructor). Hold the pad in place with tape or twist ties.

For either the simple or fractional distillation, place a boiling stone into the 10-mL round-bottom flask. Also add 8.0 mL of the unknown mixture (measured with a 10-mL graduated cylinder) to the flask. Use a hot plate and an aluminum block for heating.

Distillation

These instructions apply to both the simple and fractional distillations. Start circulating the cooling water in the condenser and adjust the heat so that the liquid boils rapidly. During the initial stages of the distillation, continue to maintain a rapid boiling rate. As the hot vapors rise, they will gradually heat up the glassware and, in the case of the fractional distillation, the fractionating column as well. Because the mass of glass and other materials is fairly large, it will take 10–20 minutes of heating before the distillation temperature begins to rise rapidly and approaches the boiling point of the distillate. (Note that this may take longer for the fractional distillation.) When the temperature begins to level off, you should soon see drops of distillate falling into the graduated cylinder.

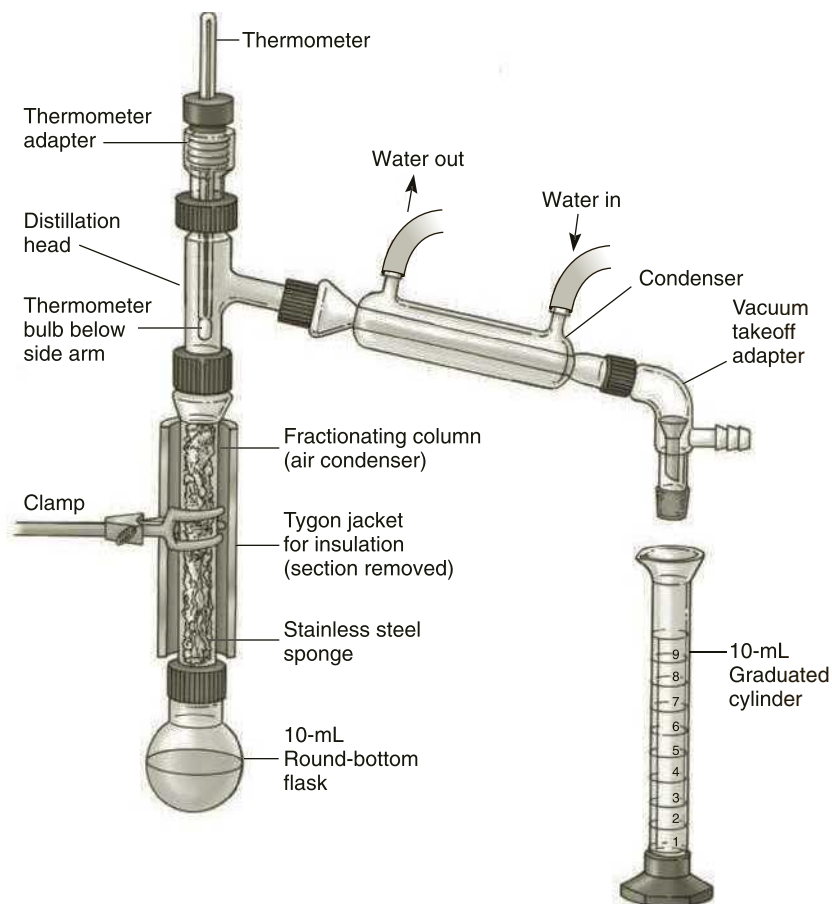

Apparatus for fractional distillation. A digital thermometer can also be used; see Section 13.4 and Figure 14.12 in the Techniques.

NOTE: For the remainder of the distillation, it is very important to regulate the temperature of the hot plate so that the distillation occurs at a rate of 1 drop per 5 seconds. **If the distillation is performed more rapidly than this, you may not achieve good separation between the liquids.** On the other hand, if the distillation is performed less rapidly than the suggested rate, the distillation temperature may be lower than it should be.

Now you will probably need to turn down the heat control to achieve the desired rate of distillation. In addition, it may be helpful to raise the round-bottom flask slightly above the aluminum block for a minute or so to cool the mixture more quickly. You should also begin recording the distillation temperature as a function of the total volume of distillate collected. Beginning at a volume of 0.5 mL, record the temperature at every 0.5-mL interval, as determined by the volume of distillate in the 10-mL graduated cylinder. After you have collected 1.0 mL of distillate, remove the 10-mL graduated cylinder and collect the next 3–4 drops of distillate in a small leak-proof vial.² Label the vial “1-mL sample.” Cap the vial tightly; otherwise, the more volatile component will evaporate more rapidly, and the composition of the mixture will change. Resume collecting the distillate in the graduated cylinder. As the distillation temperature increases, you may need to turn up the control to maintain the same rate of distillation.³ Continue to record the temperature and

² We have found GC-MS vials ideal for this purpose.

³ After the first component has distilled, it is common to observe a significant temperature drop. It may be helpful to increase the heat if this happens; however, be sure you do not exceed a distillation rate of 1 drop per 5 seconds.

volume data. When you have collected a total of 5.0 mL of distillate, take another small sample of distillate in a second small vial. (If the total volume of distillate that you can collect is less than 5.0 mL, take the last 3–4 drops.) Cap the vial and label it “5.0-mL sample.” Then continue the distillation until there is a small amount (about 0.5 mL) of liquid remaining in the distilling flask.

NOTE: Do not distill to dryness! A dry flask may crack if it is heated too hot.

The best way to stop the distillation is to turn off the hot plate and immediately raise the entire distillation apparatus off the aluminum block.

Analysis

Distillation Curve

Using the data you collected for the distillation temperature and the total volume of distillate, construct separate graphs for the simple and fractional distillations. Plot the volume in 0.5-mL increments on the x -axis and the temperature on the y -axis. Comparing the two graphs should make clear that the fractional distillation resulted in a better separation of the two liquids. Using the graph for the fractional distillation, estimate the boiling points of the two components in your mixture by noting the two regions on the graph where the temperatures leveled off. From these approximate boiling points, try to identify the two liquids in your mixture (see table at the beginning of this experiment). Note that the observed boiling point for the first component may be somewhat higher than the actual boiling point, and the observed boiling point for the second component may be somewhat lower than the actual boiling point. The reason is that the fractionating column may not be efficient enough to completely separate all of the pairs of liquids in this experiment. Therefore, it may be easier to identify the two liquids in your mixture from the gas chromatograph, as described in the next section.

Gas Chromatography

Gas chromatography is an instrumental method that separates the components of a mixture based on their boiling points. The lower-boiling component passes through the column first, followed by the higher-boiling components. The actual length of time required for a compound to pass through the column is called the **retention time** of that compound. As each component comes off the column, it is detected, and a peak is recorded that is proportional in size to the amount of the compound that was put on the column.

Gas chromatography can be used to determine the compositions of the two samples that you collected in small vials. The instructor or a laboratory assistant may either make the sample injections or allow you to make them. In the latter case, your instructor will give you adequate instructions beforehand. A reasonable sample size is 2.5 μ L. Inject the sample into the gas chromatograph and record the gas chromatogram. Depending on how effectively the two compounds were separated by the distillation, you may see one or two peaks. The lower-boiling component has a shorter retention time than the higher-boiling one. Your instructor may provide you with the actual retention times for each compound so that you can identify each peak with more certainty.

Once the gas chromatogram has been obtained, determine the relative areas of the two peaks (Technique 22, Section 22.12). You can calculate this by triangulation, or the instrument may do this electronically. In either case, you should divide each area by a response factor to account for differences in how the detector responds to

the different compounds.⁴ Calculate the percentages of the two compounds in both samples. Compare these results for the simple and fractional distillations.

8B**EXPERIMENT 8B**

Simple and Fractional Distillation (Microscale Procedure)

PROCEDURE

This experiment can also be performed using a Hickman head, although it is not as easy to monitor the volume of distillate or the rate of distillation. To perform a simple distillation, refer to Technique 14, Figure 14.7B. For a fractional distillation, see Technique 15, Figure 15.2. For both distillations, use a 10-mL round-bottom flask and a Hickman head with a side port (Technique 14, Figure 14.4B). Attach a water-cooled condenser on top of the Hickman head. It is helpful to tilt the apparatus slightly (5–10 degrees) in the direction of the side port so that the liquid in the reservoir of the Hickman head will flow toward the side port.

In both Figure 14.7B and Figure 15.2, referred to in the preceding paragraph, a thermometer is used to monitor the temperature during the distillation. It may be possible to use a stainless steel temperature probe in place of the thermometer if the temperature probe is long enough. For more discussion about the use of a temperature probe, see the first paragraph in the Apparatus section of Experiment 8A.

Follow the procedure given in Experiment 8A, except that it will be necessary to transfer the distillate from the Hickman head to a 10-mL graduated cylinder to collect data for the distillation temperature and total volume of distillate. This must be done frequently so that data can be taken at 0.5-mL intervals, as indicated in the procedure. Because you will not be able to count drops, you should try to distill at a rate of three to four minutes per mL distillate. It is important to remove as much distillate as possible each time you make a transfer. Otherwise, the next sample of distillate will be contaminated by the leftover liquid.

REPORT

Distillation Curve

Record the data for the distillation temperature as a function of the volume of distillate. Construct a graph for these data (see “Analysis” in Experiment 8A). Compare the graphs for simple and fractional distillations of the same mixture. Which distillation resulted in a better separation? Explain. Report the approximate boiling points for the two compounds in your mixture and identify the compounds.

⁴ Because response factors are instrument specific, you will be given the response factors for your instrument. Typical response factors obtained on a GowMac 69-350 gas chromatograph are hexane (1.50), cyclohexane (1.80), heptane (1.63), and toluene (1.41). These response factors were determined by injecting a mixture of equal volumes of the four liquids and determining the relative peak areas.


Gas Chromatography For both the 1-mL sample and the 5-mL sample, determine the relative areas of the two peaks, unless there is only one peak. Divide the areas by the appropriate response factors and calculate the percentage composition of the two compounds in each sample. Compare these results for the simple and fractional distillations of the same mixture. Which distillation resulted in a better separation? Explain. Identify the two compounds in your mixture. At your instructor's option, turn in the gas chromatograms with your report.

Aspirin

Aspirin is one of the most popular cure-alls available today. It is a powerful **analgesic** (relieves pain), **antipyretic** (reduces fever), **anti-inflammatory** (reduces swelling), and **antiplatelet** (slows blood-clotting) drug. Although its history as a modern medicine began only a little over a century ago, its medicinal origins actually lie in folk remedies, some of which were recognized as early as 3000 BC. Early Greek, Roman, Egyptian, Babylonian, and Chinese medical treatises recognized the ability of extracts of the willow and other salicylate-containing plants, such as meadowsweet and myrtle, to alleviate fever, pain, and inflammation. The use of meadowsweet extracts was common throughout the Middle Ages. Aspirin first appeared as a commercially available tablet in 1899. By the late 1950s, over 15 billion tablets were consumed each year. The commercial introduction of acetaminophen (Tylenol) in 1956 and of ibuprofen in 1961 caused a temporary decline in the use of aspirin. However, new uses have been found for the drug in treating heart disease (“baby aspirin”), and its popularity remains strong. Since it was first made available to the general public, it is estimated that over a trillion aspirin tablets have been consumed by patients seeking relief.


The modern history of aspirin began on June 2, 1763, when Edward Stone, a clergyman, read a paper to the Royal Society of London entitled, “An Account of the Success of the Bark of the Willow in the Cure of Agues.” By *ague*, Stone was referring to what we now call malaria, but his use of the word *cure* was optimistic; what his extract of willow bark actually did was to dramatically reduce the feverish symptoms of the disease. He was promoting his new malaria cure as a substitute for “Peruvian Bark,” an imported and expensive remedy, which we now know contains the drug quinine. Almost a century later, a Scottish physician found that Stone’s extract could also relieve the symptoms of acute rheumatism.

Soon thereafter, organic chemists working with willow bark extract and flowers of the meadowsweet plant (which gave a similar compound) isolated and identified the active ingredient as salicylic acid (from *salix*, the Latin name for the willow tree). The substance could then be chemically produced in large quantities for medical use. It soon became apparent that using salicylic acid as a remedy was severely limited by its acidic properties. The substance irritated the mucous membranes lining the mouth, esophagus, and stomach. The first attempts to circumvent this problem by using the less acidic sodium salt (sodium salicylate) were only partially successful. This substance was less irritating but had such an objectionable sweetish taste that most people could not be induced to take it. The breakthrough came at the turn of the century (1893) when Felix Hofmann, a young chemist working for the German company Bayer, devised a practical route for synthesizing acetylsalicylic acid, which was found to have all the same medicinal properties without the highly objectionable taste or the high degree of mucosal-membrane irritation. Bayer called its new product “aspirin,” a name derived from *a-* for acetyl, and the root *-spir*, from the Latin name for the meadowsweet plant, *spirea*.

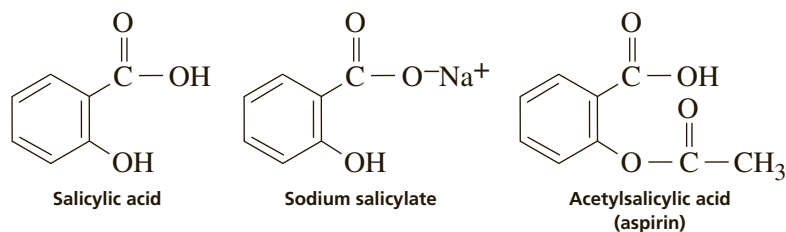


The history of aspirin is typical of many of the medicinal substances in current use. Many began as crude plant extracts or folk remedies, the active ingredients of which were isolated and their structure determined by chemists, who then improved on the original.

Through the research of J.R. Vane and others in the 1970s, aspirin's mode of action has largely been explained. A whole new class of compounds, called **prostaglandins**, has been found to be involved in the body's immune responses. Their synthesis is provoked by interference with the body's normal functioning by foreign substances or unaccustomed stimuli.

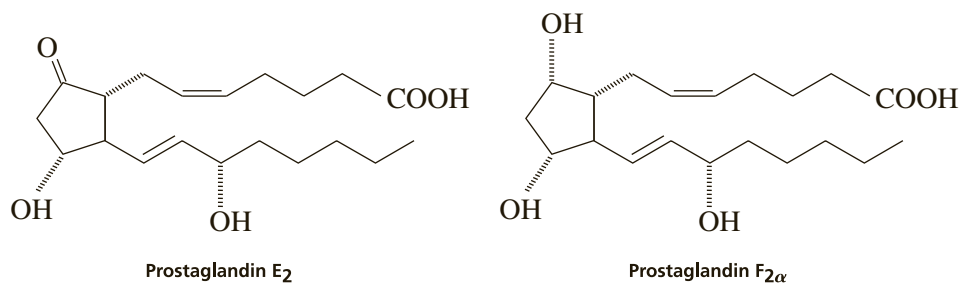

These substances are involved in a wide variety of physiological processes and are thought to be responsible for evoking pain, fever, and local inflammation. Aspirin has recently been shown to prevent bodily synthesis of prostaglandins and thus to alleviate the symptomatic portion (fever, pain, inflammation, menstrual cramps) of the body's immune responses (that is, the ones that let you know something is wrong). Research suggests that aspirin may inactivate one of the enzymes responsible for the synthesis of prostaglandins. The natural precursor for prostaglandin synthesis is **arachidonic acid**. This substance is converted to a peroxide intermediate by an enzyme called **cyclo-oxygenase**, or prostaglandin synthase. This intermediate is converted further to prostaglandin. The apparent role of aspirin is to attach an acetyl group to the active site of cyclo-oxygenase, thus rendering it unable to convert arachidonic acid to the peroxide intermediate. In this way, prostaglandin synthesis is blocked.

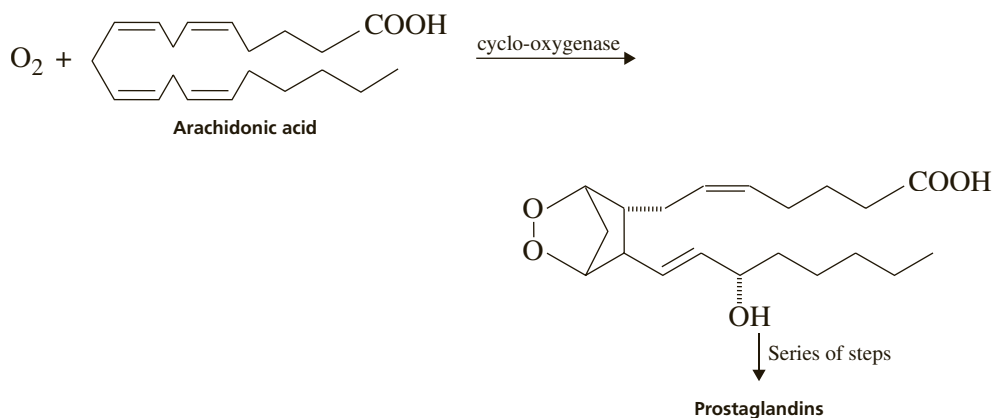

Aspirin tablets (5-grain size) are usually compounded of about 0.32 g of acetylsalicylic acid pressed together with a small amount of starch, which binds the ingredients. Buffered aspirin usually contains a basic buffering agent to reduce the acidic irritation of mucous membranes in the stomach, because the acetylated product is not totally free of this irritating effect. Bufferin contains 0.325 g of aspirin together with calcium carbonate, magnesium oxide, and magnesium carbonate as buffering agents. Combination pain relievers usually contain aspirin, acetaminophen, and caffeine. Extra- Strength Excedrin, for instance, contains 0.250 g aspirin, 0.250 g acetaminophen, and 0.065 g caffeine.

In the late 1980s scientists discovered that small daily doses of aspirin were effective in reducing the risk of blood-clotting diseases. “Baby aspirin” tablets contain about 25% (0.082 g) of the amount of acetylsalicylic acid that is contained in a regular aspirin tablet. These small tablets are often prescribed to survivors of heart attacks and strokes to prevent a recurrence. As an antiplatelet drug, aspirin prevents tiny red blood cells (platelets) from clumping together or clotting. Clotting in arteries can initiate the events that lead to arteriosclerosis. If blood clots block arteries or break loose and travel to the heart or the brain, heart attacks and strokes can occur.

Some persons are allergic to aspirin and cannot tolerate it or other salicylate-based medicines. In other people, aspirin may cause gastric irritation or ulcers and bleeding in the stomach. For this reason, doctors often prefer to prescribe acetaminophen (Tylenol). When treating children, aspirin should also be avoided in favor of Tylenol, due to a known link between aspirin consumption and Reye’s Syndrome, a disease which can be fatal; however, acetaminophen does not have any antiplatelet activity and cannot prevent or deter clotting diseases in susceptible adults. Finally, with some diseases, aspirin simply provides superior relief of pain and inflammation and is preferred over any of the newer analgesics. Following its decline in the mid-twentieth century, aspirin has undergone a resurgence and is once again a top seller in the analgesic marketplace.

REFERENCES

- Aspirin Cuts Deaths after Heart Attacks. *New Sci.* **1988**, 188 (Apr 7), 22.
- Collier, H. O. J. Aspirin. *Sci. Am.* **1963**, 209 (Nov), 96.
- Collier, H. O. J. Prostaglandins and Aspirin. *Nature* **1971**, 232 (July 2), 17.
- Disla, E.; Rhim, H. R.; Reddy, A.; and Taranta, A. Aspirin on Trial as HIV Treatment. *Nature* **1933**, 366 (Nov 18), 198.
- Jeffreys, D. *Aspirin: The Remarkable Story of a Wonder Drug*; Bloomsbury Publishing: New York, 2005.
- Kingman, S. Will an Aspirin a Day Keep the Doctor Away? *New Sci.* **1988**, 117 (Feb), 26.
- Kolata, G. Study of Reye’s-Aspirin Link Raises Concerns. *Science* **1985**, 227 (Jan 25), 391.
- Macilwain, C. Aspirin on Trial as HIV Treatment. *Nature* **1993**, 364 (Jul 29), 369.
- Nelson, N. A.; Kelly, R. C.; and Johnson, R. A. Prostaglandins and the Arachidonic Acid Cascade. *Chem. Eng. News* **1982**, 60 (Aug 16), 30.
- Pike, J. E. Prostaglandins. *Sci. Am.* **1971**, 225 (Nov), 84.
- Roth, G. J.; Stanford, N.; and Majerus, P. W. Acetylation of Prostaglandin Synthase by Aspirin. *Proc. Natl. Acad. Sci. USA* **1975**, 72, 3073.
- Street, K. W. Method Development for Analysis of Aspirin Tablets. *J. Chem. Educ.* **1988**, 65 (Oct), 914.
- Vane, J. R. Inhibition of Prostaglandin Synthesis as a Mechanism of Action for Aspirin-Like Drugs. *Nat. New Biol.* **1971**, 231 (Jun 23), 232.
- Weissmann, G. Aspirin. *Sci. Am.* **1991**, 264 (Jan), 84.

Acetylsalicylic Acid


Crystallization

Vacuum filtration

Melting point

Esterification

Aspirin (acetylsalicylic acid) can be prepared by the reaction between salicylic acid and acetic anhydride:

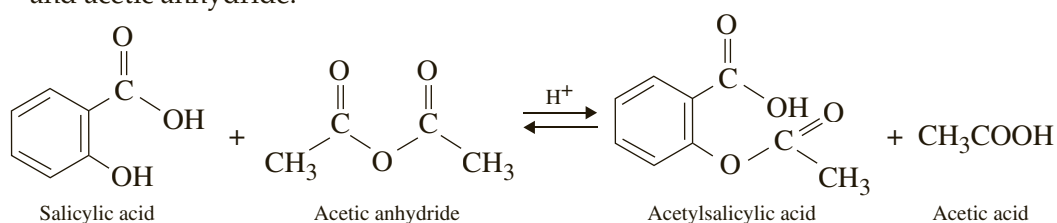

In this reaction, the **hydroxyl group** ($-\text{OH}$) on the benzene ring in salicylic acid reacts with acetic anhydride to form an **ester** functional group. Thus, the formation of acetylsalicylic acid is referred to as an **esterification** reaction. This reaction requires the presence of an acid catalyst, indicated by the H^+ above the equilibrium arrows.

When the reaction is complete, some unreacted salicylic acid and acetic anhydride will be present, along with acetylsalicylic acid, acetic acid, and the catalyst. The technique used to purify the acetylsalicylic acid from the other substances is called **crystallization**. This technique, which was introduced in Experiment 3, will be studied in more detail in Experiment 11. The basic principle is quite simple. At the end of this reaction, the reaction mixture will be hot, and all substances will be in solution. As the solution is allowed to cool, the solubility of acetylsalicylic acid will decrease, and it will gradually come out of solution, or crystallize. Because the other substances are either liquids at room temperature or are present in much smaller amounts, the crystals formed will be composed mainly of acetylsalicylic acid. Thus, a separation of acetylsalicylic acid from the other materials will have been accomplished. The purification process is facilitated by the addition of water after the crystals have formed. The water decreases the solubility of acetylsalicylic acid and dissolves some of the impurities.

The most likely impurity in the final product is salicylic acid itself, which can arise from incomplete reaction of the starting materials or from **hydrolysis** (reaction with water) of the product during the isolation steps. The hydrolysis reaction of acetylsalicylic acid produces salicylic acid. Salicylic acid and other compounds that contain a hydroxyl group on the benzene ring are referred to as **phenols**. Phenols form a highly colored complex with ferric chloride (Fe^{3+} ion). Aspirin is not a phenol, because it does not possess a hydroxyl group directly attached to the ring.

Because aspirin will not give the color reaction with ferric chloride, the presence of salicylic acid in the final product is easily detected. The purity of your product will also be determined by obtaining the melting point.

REQUIRED READING

<i>Review:</i>	Introduction to Microscale Laboratory (Experiment 1)
	Technique 8 Filtration, Sections 8.1–8.6
	Technique 9 Physical Constants, Melting Points
<i>New:</i>	Technique 5 Measurement of Volume and Weight
	Technique 6 Heating and Cooling Methods
	Technique 7 Reaction Methods, Sections 7.1–7.4
	Essay Aspirin

SPECIAL INSTRUCTIONS

This experiment involves concentrated phosphoric acid, which is highly corrosive. It will cause burns if it is spilled on the skin. Exercise care in handling it. The acetylsalicylic acid crystals should be allowed to air-dry overnight after filtration on the Hirsch funnel.

SUGGESTED WASTE DISPOSAL

Dispose of the aqueous filtrate in the container for aqueous waste.

PROCEDURE

Preparation of Acetylsalicylic Acid (Aspirin)

Prepare a hot water bath using a 250-mL beaker and a hot plate. Use about 100 mL of water and adjust the temperature to about 50°C. Weigh 0.210 g of salicylic acid ($MW = 138.1$) and place this in a **dry** 5-mL conical vial. It is not necessary for you to weigh exactly 0.210 g of salicylic acid. Try to obtain a weight within about 0.005 g of the indicated weight without spending excessive time at the balance. Record the actual weight in your notebook, and use this weight in any subsequent calculations. Using an automatic pipette or a dispensing pump, add 0.480 mL of acetic anhydride ($MW = 102.1$, $d = 1.08$ g/mL), followed by exactly one drop of concentrated phosphoric acid from a Pasteur pipette.

CAUTION



Concentrated phosphoric acid is highly corrosive. You must handle it with great care.

Add a magnetic spin vane (Technique 7, Figure 7.8A) and attach an air condenser to the vial. Clamp this assembly so that the vial is partially submerged in the hot water bath (Technique 6, Figure 6.6). Stir the mixture with the spin vane until the salicylic acid dissolves. (If the spin vane becomes stuck in the solid salicylic acid,

insert a microspatula through the air condenser into the conical vial and gently push the spin vane until it begins spinning.) Heat the mixture for 8–10 minutes after the solid dissolves to complete the reaction.

Crystallization of Acetylsalicylic Acid

Remove the vial from the water bath and allow it to cool. After the vial has cooled enough for you to handle it, detach the air condenser and remove the spin vane with forceps or a magnetic stirring bar. (If you use forceps, be sure to clean them.) Place the conical vial in a small beaker and allow the vial to cool to room temperature, during which time the acetylsalicylic acid should begin to crystallize from the reaction mixture. If it does not crystallize, scratch the walls of the vial with a glass rod (not fire-polished) and cool the mixture slightly in an ice-water bath (Technique 11, Section 11.3C) until crystallization has occurred. (Scratching the inside walls of the container often helps to initiate crystallization.) After crystal formation is complete (usually when the product appears as a solid mass), add 3.0 mL of water (measured with a 10-mL graduated cylinder) and stir thoroughly with a microspatula.

Vacuum Filtration

Set up a Hirsch funnel for vacuum filtration (see Technique 8, Section 8.3, and Figure 8.5). Moisten the filter paper with a few drops of water and turn on the vacuum (or aspirator) to the fullest extent. Transfer the mixture in the conical vial to the Hirsch funnel. When you have removed as much product as possible from the vial, add about 1 mL of cold water to the vial using a calibrated Pasteur pipette. Stir the mixture and transfer the remaining crystals and water to the Hirsch funnel. When all the crystals have been collected in the funnel, rinse them with several 0.5-mL portions of cold water. Continue drawing air through the crystals on the Hirsch funnel by suction until the crystals are nearly dry (5–10 minutes). Remove the crystals for air-drying on a watch glass or clay plate. It is convenient to hold the filter paper disc with forceps while *gently* scraping the crystals off the filter paper with a microspatula. If the paper is scraped too hard, small pieces of paper will be removed along with the crystals. To dry the crystals completely, you must set the crystals aside overnight. Weigh the dry product and calculate the percentage yield of acetylsalicylic acid ($MW = 180.2$).

Ferric Chloride Test

Test for Purity

You can perform this test on a sample of your product that is not completely dry. To determine if there is any salicylic acid remaining in your product, carry out the following procedure. Obtain three small test tubes. Add 0.5 mL of water to each test tube. Dissolve a small amount of salicylic acid in the first tube. Add a similar amount of your product to the second tube. The third test tube, which contains only solvent, will serve as the control. Add one drop of 1% ferric chloride solution to each tube and note the color after shaking. Formation of an iron–phenol complex with Fe(III) gives a definite color ranging from red to violet, depending on the particular phenol present.

Melting Point

As an additional test for purity, determine the melting point of your product (see Technique 9, Sections 9.5–9.8). The melting point must be obtained with a completely dried sample. Pure aspirin has a melting point of 135–136°C.

Place your product in a small vial, label it properly (Technique 2, Section 2.4), and submit it to your instructor.

Aspirin Tablets


Aspirin tablets are acetylsalicylic acid pressed together with a small amount of inert binding material. Common binding substances include starch, methylcellulose, and microcrystalline cellulose. You can test for the presence of starch by boiling approximately one-fourth of an aspirin tablet with 2 mL of water. Cool the liquid and add a drop of iodine solution. If starch is present, it will form a complex with the iodine. The starch–iodine complex is deep blue-violet. Repeat this test with a commercial aspirin tablet and with the acetylsalicylic acid prepared in this experiment.

QUESTIONS

1. What is the purpose of the concentrated phosphoric acid used in the first step?
2. What would happen if the phosphoric acid were left out?
3. If you used 250 mg of salicylic acid and excess acetic anhydride in the preceding synthesis of aspirin, what would be the theoretical yield of acetylsalicylic acid in moles? In milligrams?
4. What is the equation for the decomposition reaction that can occur with aspirin in water?
5. Most aspirin tablets contain five grains of acetylsalicylic acid. How many milligrams is this? (*Hint: See the essay "Aspirin."*)
6. A student performed the reaction in this experiment using a water bath at 90°C instead of 50°C. The final product was tested for the presence of phenols with ferric chloride. This test was negative (no color observed); however, the melting point of the dry product was 122–125°C. Explain these results as completely as possible.
7. If the aspirin crystals were not completely dried before the melting point was determined, what effect would this have on the observed melting point?


Analgesics

Acylated aromatic amines (those having an acyl group, $\text{R}-\overset{\text{O}}{\parallel}{\text{C}}-$, substituted on nitrogen) are important in over-the-counter headache remedies. Over-the-counter drugs are those you may buy without a prescription. Acetanilide, phenacetin, and acetaminophen are mild analgesics (relieve pain) and antipyretics (reduce fever) and are important, along with aspirin, in many nonprescription drugs.

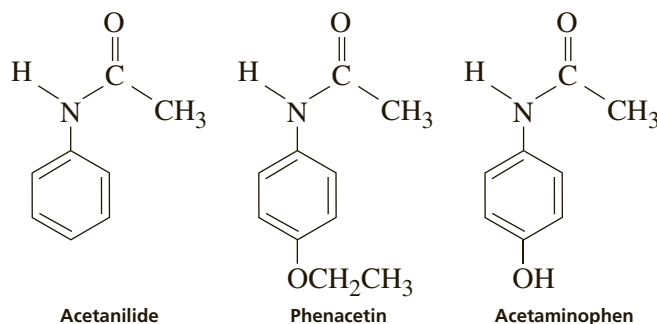


The discovery that acetanilide was an effective antipyretic came about by accident in 1886. Two doctors, Cahn and Hepp, had been testing naphthalene as a possible **vermifuge** (an agent that expels worms). Their early results on simple worm cases were very discouraging, so Dr. Hepp decided to test the compound on a patient with a larger variety of complaints, including worms—a sort of shotgun approach. A short time later, Dr. Hepp excitedly reported to his colleague, Dr. Cahn, that naphthalene had miraculous fever-reducing properties.

In trying to verify this observation, the doctors discovered that the bottle they thought contained naphthalene apparently did not. In fact, the bottle brought to them by their assistant had a label so faint as to be illegible. They were sure that the sample was not naphthalene, because it had no odor. Naphthalene has a strong odor reminiscent of mothballs. So close to an important discovery, the doctors were nevertheless stymied. They appealed to Hepp's cousin, who was a chemist in a nearby dye factory, to help them identify the unknown compound. This compound turned out to be acetanilide, a compound with a structure not at all like that of naphthalene. Certainly, Hepp's unscientific and risky approach would be frowned on by doctors today; and to be sure, the Food and Drug Administration (FDA) would never allow human testing before extensive animal testing (consumer protection has greatly progressed). Nevertheless, Cahn and Hepp made an important discovery.

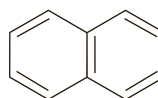


Naphthalene

In another instance of serendipity, Cahn and Hepp's publication, describing their experiments with acetanilide, caught the attention of Carl Duisberg, director of research at the Bayer Company in Germany. Duisberg was confronted with the problem of how to profitably get rid of nearly 50 tons of *p*-aminophenol, a by-product of the synthesis of one of Bayer's other commercial products. He immediately saw the possibility of converting *p*-aminophenol to a compound similar in structure to acetanilide by putting an acyl group on the nitrogen. It was then believed, however, that all compounds having a hydroxyl group on a benzene ring (that is, phenols) were toxic. Duisberg devised a scheme of structural modification of *p*-aminophenol to synthesize the compound phenacetin. The reaction scheme is shown here.

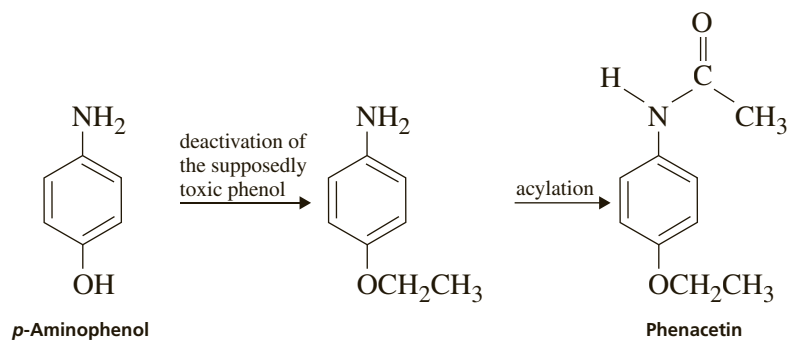


Phenacetin turned out to be a highly effective analgesic and antipyretic. A common form of combination pain reliever called an APC tablet was once available. An APC tablet contained Aspirin, Phenacetin, and Caffeine (hence, APC). Phenacetin is no longer used in commercial pain-relief preparations as it was discovered that not all aromatic hydroxyl groups lead to toxic compounds. Today the compound acetaminophen is very widely used as an analgesic in place of phenacetin.

Another analgesic, structurally similar to aspirin, that has found some application is **salicylamide**. Salicylamide is an ingredient in some pain-relief preparations, although its use is declining.

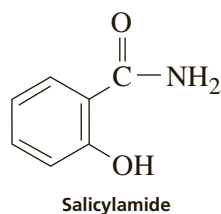

Upon continued or excessive use, acetanilide can cause a serious blood disorder called **methemoglobinemia**. In this disorder, the central iron atom in hemoglobin is converted from Fe(II) to Fe(III) to give methemoglobin. Methemoglobin will not function as an oxygen carrier in the bloodstream. The result is a type of anemia (deficiency of hemoglobin or lack of red blood cells). Phenacetin and acetaminophen cause the same disorder, but to a much lesser degree. Because they are also more effective as antipyretic and analgesic drugs than acetanilide, they are preferred remedies. Acetaminophen is marketed under a variety of trade names, including Tylenol, Datril, and Panadol, and is often successfully used by people who are allergic to aspirin.

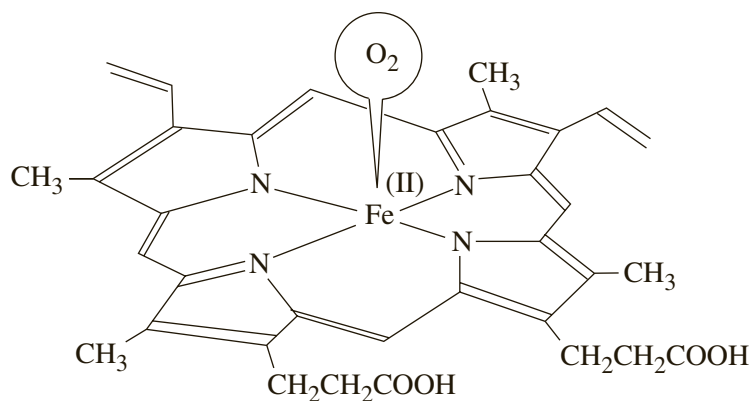

Heme portion of blood-oxygen carrier, hemoglobin.


More recently, a new drug has appeared in over-the-counter preparations. This drug is **ibuprofen**, which was initially marketed as a prescription drug in the United States under the name Motrin. Ibuprofen was first developed and patented in England in 1961. The United States obtained marketing rights in 1974. Ibuprofen is now sold without a prescription under several brand names, which include Advil, Motrin, and Nuprin. Ibuprofen is principally an anti-inflammatory drug, but it is also effective as an analgesic and an antipyretic. It is particularly effective in treating the symptoms of rheumatoid arthritis and menstrual cramps. Ibuprofen appears to control the production of prostaglandins, which parallels aspirin's mode of action. An important advantage of ibuprofen is that it is a very powerful pain reliever. One 200-mg tablet is as effective as two tablets (650 mg) of aspirin. Furthermore, ibuprofen has a more advantageous dose-response curve, which means that taking two tablets of this drug is approximately twice as effective as one tablet for certain types of pain. Aspirin and acetaminophen reach their maximum effective dose at two tablets. Little additional relief is gained at doses above that level. Ibuprofen, however, continues to increase its effectiveness up to the 400-mg level (the equivalent of four tablets of aspirin or acetaminophen). Ibuprofen is a relatively safe drug, but its use should be avoided in cases of aspirin allergy, kidney problems, ulcers, asthma, hypertension, or heart disease.

Analgesics and caffeine in some common preparations

	<i>Aspirin</i>	<i>Acetaminophen</i>	<i>Caffeine</i>
Aspirin*	0.325 g	—	—
Anacin	0.400 g	—	0.032 g
Bufferin	0.325 g	—	—
Cope	0.421 g	—	0.032 g
Excedrin (Extra-Strength)	0.250 g	0.250 g	0.065 g
Tylenol	—	0.325 g	—
B. C. Tablets	0.325 g	—	0.016 g
Advil	—	—	—
Aleve	—	—	—
Orudis	—	—	—


Note: Nonanalgesic ingredients (e.g., buffers) are not listed.

*5-grain tablet (1 grain = 0.0648 g).

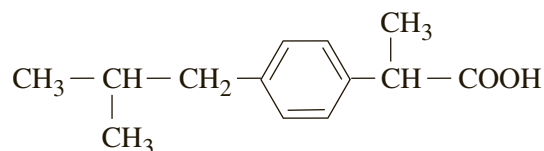

Ibuprofen

The Food and Drug Administration has also approved two other drugs with similar structures to ibuprofen for over-the-counter use as pain relievers. These new drugs are known by their generic names, **naproxen** and **ketoprofen**. Naproxen is often administered in the form of its sodium salt. Naproxen and ketoprofen can be used to alleviate the pain of headaches, toothaches, muscle aches, backaches, arthritis, and menstrual cramps, and they can also be used to reduce fever. They appear to have a longer duration of action than the older analgesics.

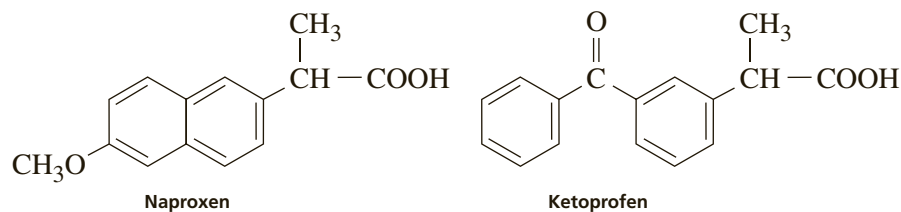

Salicylamide	Ibuprofen	Ketoprofen	Naproxen
—	—	—	—
—	—	—	—
—	—	—	—
—	—	—	—
—	—	—	—
0.095 g	—	—	—
—	0.200 g	—	—
—	—	—	0.220 g
—	—	0.0125 g	—

REFERENCES

- Bugg, C. E.; Carson, W. M.; and Montgomery, J. A. Drugs by Design. *Sci. Am.* **1993**, 269 (Dec), 92.
- Grosser, T.; Smyth, E. M.; and Garret, A. F. Anti-Inflammatory, Antipyretic and Analgesic Agents. In *Goodman and Gilman's: The Pharmacological Basis of Therapeutics*, 12th ed.; Bronton, L., Chabner, B., and Knollman, B., Eds.; Macmillan: New York, 2010.
- Hansch, C. Drug Research or the Luck of the Draw. *J. Chem. Educ.* **1974**, 51, 360.
- Hart, C.; and Ksir, C. Internal Analgesics. *Drugs, Society, and Human Behavior*, 15th ed.; McGraw-Hill: New York, 2012.
- The New Pain Relievers. *Consum. Rep.* **1984**, 49 (Nov), 636–638.
- Senozan, N. M. Methemoglobinemia: An Illness Caused by the Ferric State. *J. Chem. Educ.* **1985**, 62 (Mar), 181.
- Tietze, K. L., Ed. Section II, "Pain and Fever Disorders." In *Handbook of Non-Prescription Drugs*, 17th ed.; Krinsky, D. L., Ed.; American Pharmacists Association: Washington, DC, 2012.


Isolation of the Active Ingredient in an Analgesic Drug

Extraction


Filtration

Melting point


Most analgesic (pain-relieving) drugs found on the shelves of any drug or grocery store generally fall into one of four categories. These drugs may contain **acetylsalicylic acid**, **acetaminophen**, or **ibuprofen** as the active ingredient, or some **combination** of these compounds may be used in a single preparation. All tablets, regardless of type, contain a large amount of starch or other inert substance. This material acts as a binder to keep the tablet from falling apart and to make it large enough to handle. Some analgesic drugs also contain caffeine or buffering agents. In addition, many tablets are coated to make them easier to swallow and to prevent users from experiencing the unpleasant taste of the drugs.

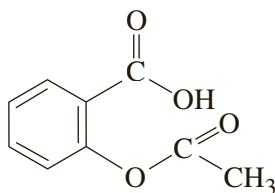

Acetylsalicylic acid

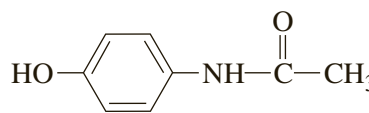

Acetaminophen

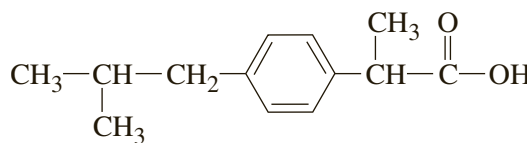

Ibuprofen

The three drugs, along with their melting points (MP) and common brand names, follow:

<i>Drug</i>	<i>MP</i>	<i>Brand Names</i>
Acetylsalicylic acid	135–136°C	Aspirin, ASA, acetylsalicylic acid, generic aspirin, Empirin
Acetaminophen	169–170.5°C	Tylenol, Datril, Panadol, nonaspirin pain reliever (various brands)
Ibuprofen	75–77°C	Advil, Brufen, Motrin, Nuprin

The purpose of this experiment is to demonstrate some important techniques that are applied throughout this textbook and to allow you to become accustomed to working in the laboratory at the microscale level. More specifically, you will extract

(dissolve) the active ingredient of an analgesic drug by mixing the powdered tablet with a solvent, methanol. Two steps are required to remove the fine particles of binder, which remain suspended in the solvent. First, you will use centrifugation to remove most of the binder. The second step will be a filtration technique using a Pasteur pipette packed with alumina (finely ground aluminum oxide). The solvent will then be evaporated to yield the solid analgesic, which will be collected by filtration on a Hirsch funnel. Finally, you will test the purity of the drug by doing a melting-point determination.

REQUIRED READING

<i>Review:</i>	Experiment 1	Introduction to Microscale Laboratory
<i>New:</i>	Technique 7	Reaction Methods, Section 7.9
	Technique 8	Filtration, Sections 8.1–8.6
	Technique 9	Physical Constants of Solids: The Melting Point

SPECIAL INSTRUCTIONS

You will be allowed to select an analgesic that is a member of one of the categories described previously. You should use an uncoated tablet that contains only a single ingredient analgesic and binder. If it is necessary to use a coated tablet, try to remove the coating when the tablet is crushed. To avoid decomposition of aspirin, it is essential to minimize the length of time that it remains dissolved in methanol. Do not stop this experiment until after the drug is dried on the Hirsch funnel.

SUGGESTED WASTE DISPOSAL

Dispose of any remaining methanol in the waste container for nonhalogenated organic solvents. Place the alumina in the container designated for wet alumina.


PROCEDURE

Extraction of Active Ingredient

If you are isolating aspirin or acetaminophen, use *one* tablet in this procedure. If you are isolating ibuprofen, use *two* tablets. Using a pestle, crush the tablet (or tablets) between two pieces of weighing paper. If the tablet is coated, try to remove fragments of the coating material with forceps after the tablet is first crushed. Add all the powdered material to a 3-mL conical vial. Using a calibrated Pasteur pipette, add about 2 mL of methanol to the vial. Cap the vial and mix thoroughly by shaking. Loosen the cap at least once during the mixing process to release any pressure that may build up in the vial.

Allow the undissolved portion of the powder to settle in the vial. A cloudy suspension may remain even after 5 minutes or more. You should wait only until it is obvious that the larger particles have settled completely. Using a filter-tip pipette (Technique 8, Figure 8.9), transfer the liquid phase to a centrifuge tube. Add a second 2-mL portion of methanol to the conical vial and repeat the shaking process described previously. After the solid has settled, transfer the liquid phase to the centrifuge tube containing the first extract.

Place the tube in a centrifuge along with another centrifuge tube of equal weight on the opposite side. Centrifuge the mixture for two to three minutes. The suspended solids should collect on the bottom of the tube, leaving a clear or nearly clear **supernatant liquid**, the liquid above the solid. If the liquid is still quite cloudy, repeat the centrifugation for a longer period or at a higher speed. Being careful not to disturb the solid at the bottom of the tube, transfer the supernatant liquid with a Pasteur pipette to a test tube or small beaker.

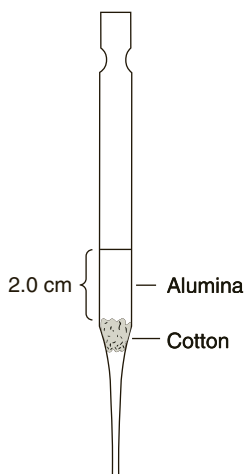

Column for purifying an analgesic drug.

Column Chromatography

Prepare an alumina column using a Pasteur pipette, as shown in the figure. Insert a small ball of cotton into the top of the column. Using a long, thin object such as a glass stirring rod or a wooden applicator stick, push the cotton down so that it fits into the Pasteur pipette where the constriction begins. Add about 0.5 g of alumina to the pipette and tap the column with your finger to pack the alumina. Clamp the pipette in a vertical position so that the liquid can drain from the column into a small beaker or a 5-mL conical vial. Place a small beaker under the column. With a calibrated Pasteur pipette, add about 2 mL of methanol to the column and allow the liquid to drain until the level of the methanol just reaches the top of the alumina. Once methanol has been added to the alumina, the top of the alumina in the column should not be allowed to run dry. If necessary, add more methanol.

NOTE: It is essential that the methanol not be allowed to drain below the surface of the alumina.

When the level of the methanol reaches the surface of the alumina, transfer the solution containing the drug from the beaker or test tube to the column using a Pasteur pipette. Collect the liquid that passes through the column into a 5-mL conical vial. When all the liquid from the beaker has been added to the column and has penetrated the alumina, add an additional 1 mL of methanol to the column and allow it to drain. This ensures that all the analgesic drug has been eluted from the column.

Evaporation of Solvent

If you are isolating aspirin, it is essential that the following evaporation procedure be completed in 10–15 minutes. Otherwise, the aspirin may partially decompose. Using a Pasteur pipette, transfer about half the liquid in the 5-mL conical vial to another small container. Evaporate the methanol in the 5-mL conical vial using a water bath at about 50°C (Technique 7, Section 7.10).¹ To speed evaporation, direct a gentle stream of dry air or nitrogen into the vial containing the liquid. Evaporate the solvent until the volume is less than about 1 mL. Then add the remainder of the liquid and continue evaporation.

When the solvent has completely evaporated or it is apparent that the remaining liquid is no longer evaporating, remove the vial from the water bath (or sand bath) and allow it to cool to room temperature. (The volume of liquid should be less than 0.5 mL when you discontinue evaporation.) If liquid remains, which is likely with the ibuprofen- or acetaminophen-containing analgesics, place the cool vial in an ice-water bath for 10–15 minutes. Prepare the ice-water bath in a small beaker, using both ice and water. Be sure that the vial cannot tip over. Crystallization of the product may occur more readily if you scrape the inside of the vial with a microspatula or a glass rod (not fire-polished). If the solid is hard and clumped,

¹ As an alternative, you may use a sand bath at about 50°C.

you should use a microspatula to break up the solid as much as possible before going on to the next step.

Vacuum Filtration

Set up a Hirsch funnel for vacuum filtration (see Technique 8, Section 8.3, and Figure 8.5). Moisten the filter paper with a few drops of methanol and turn on the vacuum (or aspirator) to the fullest extent. Use a microspatula to transfer the material in the conical vial to the Hirsch funnel. The vacuum will draw any remaining solvent from the crystals. Allow the crystals to dry for 5–10 minutes while air is drawn through the crystals in the Hirsch funnel.

Carefully scrape the dried crystals from the filter paper onto a tared (previously weighed) watch glass. If necessary, use a spatula to break up any remaining large pieces of solid. Allow the crystals to air-dry on the watch glass. To determine when the crystals are dry, move them around with a dry spatula. When the crystals no longer clump or cling to the spatula, they should be dry. If you are working with ibuprofen, the solid will be slightly sticky even when it is completely dried. Weigh the watch glass with the crystals to determine the weight of analgesic drug that you have isolated. Use the weight of the active ingredient specified on the label of the container as a basis for calculating the weight percentage recovery.

Use a small sample of the crystals to determine the melting point (see Technique 9, Sections 9.5–9.8). Crush the crystals into a powder, using a stirring rod, in order to determine their melting point. You may observe some “sweating” or shrinkage (see Technique 9, Section 9.7) before the substance actually begins to melt. The beginning of the melting-point range is when actual melting is observed, not when the solid takes on a slightly wet or shiny appearance or when shrinkage occurs. If you have isolated ibuprofen, the melting point may be somewhat lower than that given at the beginning of this experiment.

At the instructor’s option, place your product in a small vial, label it properly (Technique 2, Section 2.4), and submit it to your instructor.

QUESTIONS

1. Why was the percentage recovery less than 100%? Give several reasons.
2. Why was the tablet crushed?
3. What was the purpose of the centrifugation step?
4. What was the purpose of the alumina column?
5. If 185 mg of acetaminophen were obtained from a tablet containing 350 mg of acetaminophen, what would be the weight percentage recovery?
6. A student who was isolating aspirin stopped the experiment after the filtration step with alumina. One week later, the methanol was evaporated and the experiment was completed. The melting point of the aspirin was found to be 110–115°C. Explain why the melting point was low and why the melting range was so wide.

Acetaminophen

Decolorization


Filtration

Crystallization

Use of a Craig tube or Hirsch funnel

Preparation of an amide

Preparation of acetaminophen involves treating an amine with an acid anhydride to form an amide. In this case, *p*-aminophenol, the amine, is treated with acetic anhydride to form acetaminophen (*p*-acetamidophenol), the amide.

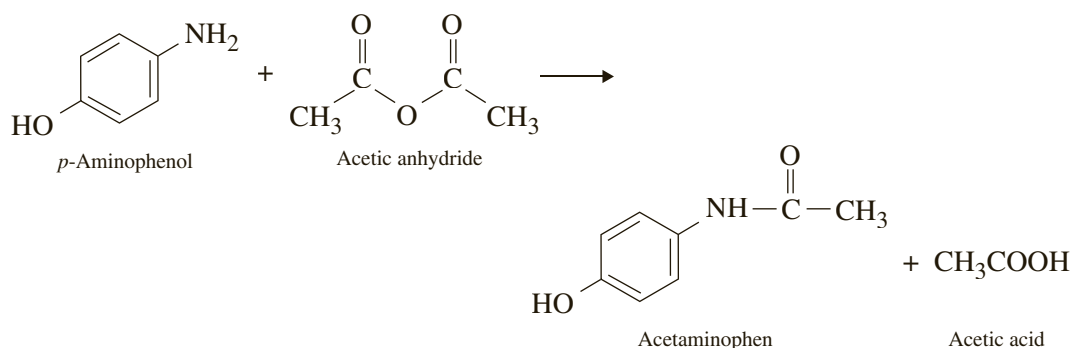

The crude solid acetaminophen contains dark impurities carried along with the *p*-aminophenol starting material. These impurities, which are dyes of unknown structure, are formed from oxidation of the starting phenol. Although the amount of the dye impurity is small, it is intense enough to impart color to the crude acetaminophen. Most of the colored impurity is destroyed by heating the crude product with sodium dithionite (sodium hydrosulfite $\text{Na}_2\text{S}_2\text{O}_4$). The dithionite reduces double bonds in the colored dye to produce colorless substances.

The decolorized acetaminophen is collected on a Hirsch funnel. It is further purified by crystallization from a methanol/water mixture. There are two procedures given in this experiment. Experiment 11A involves crystallization using a Craig tube, whereas Experiment 11B is a larger-scale reaction involving an Erlenmeyer flask and Hirsch funnel for crystallization.

REQUIRED READING

<i>Review:</i> Experiment 1	Introduction to Microscale Laboratory (Experiment 1)
Techniques 5 and 6	
Technique 7	Reaction Methods, Sections 7.1–7.3

	Technique 8	Filtration, Sections 8.1–8.6
	Technique 9	Physical Constants, Melting Points
<i>New:</i>	Technique 8	Filtration, Section 8.7
	Technique 11	Crystallization
	Essay	Analgesics

SPECIAL INSTRUCTIONS

Acetic anhydride can cause irritation of tissue, especially in nasal passages. Avoid breathing the vapor, and avoid contact with skin and eyes. *p*-Aminophenol is a skin irritant and is toxic.

WASTE DISPOSAL

Aqueous solutions obtained from filtration operations should be poured into the container designated for aqueous wastes. This includes the filtrates from the methanol and water crystallization steps.

NOTES TO THE INSTRUCTOR

The *p*-aminophenol acquires a black color on standing due to air oxidation. It is best to use a recently purchased sample, which usually is gray. If necessary, black material can be decolorized by heating it in a 10% aqueous solution of sodium dithionite (sodium hydrosulfite) prior to starting the experiment.

11A EXPERIMENT 11A

Acetaminophen (Microscale Procedure)

PROCEDURE

Reaction Mixture

Weigh about 0.150 g of *p*-aminophenol ($MW = 109.1$) and place this in a 5-mL conical vial. Using an automatic pipette (or a dispensing pump or a graduated pipette), add 0.450 mL of water and 0.165 mL of acetic anhydride ($MW = 102.1$, $d = 1.08$ g/mL). Place a spin vane in the conical vial and attach an air condenser.

Heating

Heat the reaction mixture with an aluminum block or sand bath at about 120°C (see inset in Technique 6, Figure 6.2A) and stir gently. If you are using a sand bath, the conical vial should be partially buried in the sand so that the vial is nearly at the bottom of the sand bath. After the solid has dissolved (it may dissolve, precipitate, and redissolve), heat the mixture for an additional 20 minutes to complete the reaction.

Isolation of Crude Acetaminophen

Remove the vial from the heat and allow it to cool. When the vial has cooled to the touch, detach the air condenser and remove the spin vane with clean forceps or a magnet. Rinse the spin vane with two or three drops of warm water, allowing the water to drop into the conical vial. Place the conical vial in a small beaker and let it cool to room temperature. If crystallization has not occurred, scratch the inside of the vial with a glass stirring rod to initiate crystallization. Cool the mixture thoroughly in an ice bath for 15–20 minutes and collect the crystals by vacuum filtration on a Hirsch funnel (see Technique 8, Section 8.3, and Figure 8.5). Rinse the vial with about 0.5 mL of ice water and transfer this mixture to the Hirsch funnel. Wash the crystals on the funnel with two additional 0.5-mL portions of ice water. Dry the crystals for 5–10 minutes by allowing air to be drawn through them while they remain on the Hirsch funnel. Transfer the product to a watch glass or clay plate and allow the crystals to dry in air. It may take several hours for the crystals to dry completely, but you may go on to the next step before they are totally dry. Weigh the crude product and set aside a small sample for a melting-point determination and a color comparison after the next step. Calculate the percentage yield of crude acetaminophen ($MW = 151.2$). Record the appearance of the crystals in your notebook.

Decolorization of Crude Acetaminophen

Dissolve 0.2 g of sodium dithionite (sodium hydrosulfite) in 1.5 mL of water in a 5-mL conical vial. Add your crude acetaminophen to the vial. Heat the mixture at about 100°C for 15 minutes, with occasional stirring with a microspatula. Some of the acetaminophen will dissolve during the decolorization process. Cool the mixture thoroughly in an ice bath for about 10 minutes to reprecipitate the decolorized acetaminophen (scratch the inside of the vial, if necessary, to induce crystallization). Collect the purified material by vacuum filtration on a Hirsch funnel, using small portions (about 0.5 mL total) of ice water to aid the transfer. Dry the crystals for 5–10 minutes by allowing air to be drawn through them while they remain on the Hirsch funnel. You may go on to the next step before the material is totally dry. Weigh the purified acetaminophen and compare the color of the purified material to that obtained in the preceding paragraph.

Crystallization of Acetaminophen

Place the purified acetaminophen in a Craig tube. Crystallize the material from a solvent mixture composed of 50% water and 50% methanol by volume (aluminum block or sand bath set at about 100°C). Follow the crystallization procedure described in Technique 11, Section 11.4, and Figure 11.6. The solubility of acetaminophen in this hot (nearly boiling) solvent is about 0.2 g/mL. Although you can use this as a rough indication of how much solvent is required to dissolve the solid, you should still use the technique shown in Figure 11.6 to determine how much solvent to add. Add small portions (several drops) of hot solvent until the solid dissolves. Step 2 in Figure 11.6 (removal of insoluble impurities) should not be required in this crystallization. When the solid has dissolved, place the Craig tube in a 10-mL Erlenmeyer flask, insert the inner plug of the Craig tube, and let the solution cool.

When the mixture has cooled to room temperature, place the Craig tube in an ice-water bath for several minutes. If necessary, induce crystallization by gently scratching the inside of the Craig tube with your microspatula (Technique 11, Section 11.8B). Because acetaminophen may crystallize *slowly* from the solvent,

continue to cool the Craig tube in an ice bath for at least 10 minutes. Collect the crystals using the apparatus shown in Technique 8, Figure 8.11. Place the assembly in a centrifuge (be sure it is balanced by a centrifuge tube filled with water so that both tubes contain the same weight) and turn on the centrifuge for several minutes. Collect the crystals on a watch glass or piece of smooth paper, as shown in Technique 11, Figure 11.6. Set the crystals aside to air-dry. Very little additional time should be required to complete the drying.

Yield Calculation and Melting-Point Determination

Weigh the crystallized acetaminophen ($MW = 151.2$) and calculate the percentage yield. This calculation should be based on the number of moles of the limiting reagent used at the beginning of this procedure. Determine the melting point of the product. Compare the melting point of the final product with that of the crude acetaminophen. Also compare the colors of the crude, decolorized, and pure acetaminophen. Pure acetaminophen melts at $169.5\text{--}171^\circ\text{C}$. Place your product in a properly labeled vial and submit it to your instructor.

11B EXPERIMENT 11B

Acetaminophen (Semimicroscale Procedure)

PROCEDURE

Reaction Mixture

Weigh about 0.400 g of *p*-aminophenol ($MW = 109.1$) and place this in a 5-mL conical vial. Using an automatic pipette (or a dispensing pump or a graduated pipette), add 1.20 mL of water and 0.450 mL of acetic anhydride ($MW = 102.1$, $d = 1.08$ g/mL). Place a spin vane in the conical vial and attach an air condenser.

Heating

Heat the reaction mixture with an aluminum block or a sand bath at about 120°C (see inset in Technique 7, Figure 7.6A) and stir gently. If you are using an aluminum block, position the vial so that it is touching the surface of the hot plate and place aluminum collars around the vial. If you are using a sand bath, the conical vial should be partially buried in the sand so that the vial is nearly at the bottom of the sand bath. After the solid has dissolved (it may dissolve, precipitate, and redissolve), heat the mixture for an additional 20 minutes to complete the reaction.

Isolation of the Crude Acetaminophen

Remove the vial from the heat and allow it to cool. When the vial has cooled to the touch, detach the air condenser and remove the spin vane with clean forceps or a magnet. Rinse the spin vane with two or three drops of warm water, allowing the water to drop into the conical vial. Place the conical vial in a small beaker and let it cool to room temperature. If crystallization has not occurred, scratch the inside of the vial with a glass stirring rod to initiate crystallization. Cool the mixture thoroughly in an ice bath for 15–20 minutes and collect the crystals by vacuum filtration on a Hirsch funnel (see Technique 8, Section 8.3, and Figure 8.5). Rinse the vial with about 0.5 mL of ice-cold water and transfer this mixture to the Hirsch funnel.

Repeat this rinsing with two additional 0.5-mL portions of ice-cold water. Dry the crystals for 5–10 minutes by allowing air to be drawn through them while they remain on the Hirsch funnel. Transfer the product to a watch glass or a clay plate and allow the crystals to dry in air. It may take several hours for the crystals to dry completely, but you may go on to the next step before they are totally dry. Weigh the crude product and set aside a small sample for a melting point determination and a color comparison after the next step. Calculate the percentage yield of crude acetaminophen ($MW = 151.2$). Record the appearance of the crystals in your notebook.

Decolorization of Crude Acetaminophen

Dissolve 0.5 g of sodium dithionite (sodium hydrosulfite) in 4.0 mL of water in a small Erlenmeyer flask. Add your crude acetaminophen to the flask. Heat the mixture at about 100°C for 15 minutes, with occasional swirling. Some of the acetaminophen will dissolve during the decolorization process. Cool the mixture thoroughly in an ice-water bath for about 10 minutes to reprecipitate the decolorized acetaminophen (scratch the inside of the flask, if necessary, to induce crystallization). Collect the purified material by vacuum filtration on a Hirsch funnel using small portions (about 1.0 mL total) of ice-cold water to aid the transfer. Dry the crystals for 5–10 minutes by allowing air to be drawn through them while they remain on the Hirsch funnel. You may go on to the next step before the material is totally dry. Weigh the purified acetaminophen and compare the color of the purified material to that obtained earlier.

Crystallization of Acetaminophen

Follow the semimicroscale crystallization procedure described in Technique 11, Section 11.3, and shown in Figure 11.4. Step 2 in Figure 11.4 (removal of insoluble impurities) will not be required in this crystallization.

Place all your acetaminophen in a 10-mL Erlenmeyer flask. In another flask, place about 3 mL of a solvent mixture composed of 50% water and 50% methanol, *with a boiling stone*, and put it on the hot plate. When the solvent begins to boil, start adding the hot solvent slowly to the acetaminophen using a Pasteur pipette. At this point, place both flasks on the hot plate to keep them hot. Continue to add the boiling solvent to the flask containing the acetaminophen until the solid *just dissolves*. Because the solubility of acetaminophen in this nearly boiling solvent is only about 0.2 g/mL, you will likely not use all the boiling solvent. The idea is to add the *minimum* amount of boiling solvent that *just dissolves* the solid.

Once the solid is dissolved, cork the flask and allow the contents of the flask to cool slowly to room temperature. Pure acetaminophen should crystallize out of the solvent. If solid does not form, scratch the inside of the flask with your microspatula. Place the flask in an ice bath to complete the crystallization for at least 10 minutes. Transfer the solid from the flask to a Hirsch funnel (Technique 8, Section 8.3, and Figure 8.5). Rinse the flask with about 0.5 mL of ice-cold solvent (50% methanol/50% water) and transfer this mixture to the Hirsch funnel. Repeat this rinsing with an additional 0.5-mL portion of ice-cold solvent. Dry the crystals for 5–10 minutes by allowing air to be drawn through them while they remain on the Hirsch funnel. Transfer the product to a watch glass or a clay plate and allow the crystals to dry in air. Let the crystals dry until the next laboratory period.


Yield Calculation and Melting Point Determination

Weigh the crystallized acetaminophen and calculate the percentage yield ($MW = 151.2$). This calculation should be based on the number of moles of the limiting reagent used at the beginning of this procedure. Determine the melting point of the product.

Compare the melting point of the final product with that of the crude acetaminophen. Also compare the colors of the crude, decolorized, and pure acetaminophen. Pure acetaminophen melts at 169.5–171°C. Place your product in a properly labeled vial and submit it to your instructor.

QUESTIONS

1. During the crystallization of acetaminophen, why was the mixture cooled in an ice bath?
2. In the reaction between *p*-aminophenol and acetic anhydride to form acetaminophen, 0.450 mL of water was added. What was the purpose of the water?
3. Why should you use a minimum amount of water to rinse the conical vial while transferring the purified acetaminophen to the Hirsch funnel?
4. If 0.130 g of *p*-aminophenol is allowed to react with excess acetic anhydride, what is the theoretical yield of acetaminophen in moles? In grams?
5. Give two reasons, discussed in Experiments 9 and 11, why the crude product in most reactions is not pure.
6. Phenacetin has the structure shown. Write an equation for its preparation, starting from 4-ethoxyaniline.

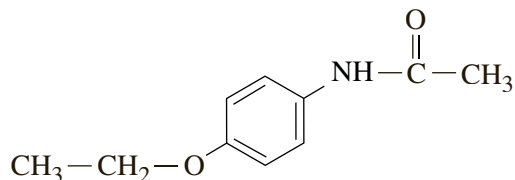

Identification of Drugs

Frequently, a chemist is called on to identify a particular unknown substance. If there is no prior information to work from, this can be a formidable task. There are several million known compounds, both inorganic and organic. For a completely unknown substance, the chemist must often use every available method. If the unknown substance is a mixture, then the mixture must be separated into its components and each component identified separately. A pure compound can often be identified from its physical properties (melting point, boiling point, density, refractive index, and so on) and a knowledge of its functional groups. These groups can be identified by the reactions that the compound is observed to undergo or by spectroscopy (infrared, ultraviolet, nuclear magnetic resonance, and mass spectroscopy). The techniques necessary for this type of identification are introduced in a later section.

A somewhat simpler situation often arises in drug identification. The scope of drug identification is more limited, and the chemist working in a hospital trying to identify the drug in an overdose or the law enforcement officer trying to identify a suspected illicit drug or a poison usually has some prior clues to work from. So does the medicinal chemist working for a pharmaceutical manufacturer who might be trying to discover why a competitor's product may be better.

Consider a drug overdose case as an example. The patient is brought into the emergency ward of a hospital. This person may be in a coma or a hyperexcited state, have an allergic rash, or clearly be hallucinating. These physiological symptoms are themselves a clue to the nature of the drug. Samples of the drug may be found in the patient's possession. Correct medical treatment may require a rapid and accurate identification of a drug powder or capsule. If the patient is conscious, the necessary information can be elicited orally; if not, the drug must be examined. If the drug is in the form of a tablet or capsule, the process is often simple because many drugs are coded by a manufacturer's trademark or logo, by shape (round, oval, or bullet shape), by formulation (tablet, gelatin capsule, or time-release microcapsule), and by color. Some drugs also bear an imprinted number or code.

It is more difficult to identify a powder, but such identification may be easy under some circumstances. Plant drugs are often easily identified because they contain microscopic bits and pieces of the plant from which they are obtained. This cellular debris is often characteristic for certain types of drugs, and they can be identified on this basis alone. A microscope is all that is needed. Sometimes chemical color tests can be used as confirmation. Certain drugs give rise to characteristic colors when treated with special reagents. Other drugs form crystalline precipitates of characteristic color and crystal structure when treated with appropriate reagents.

If the drug itself is not available and the patient is unconscious (or dead), identification may be more difficult. It may be necessary to pump the stomach or bladder contents of the patient (or corpse) or to obtain a blood sample. These samples of stomach fluid, urine, or blood would be extracted with an appropriate organic solvent, and the extract would be analyzed.

Often the final identification of a drug, as extracted from stomach fluid, urine, or blood hinges on some type of **chromatography**. Thin-layer chromatography (TLC) is often used. Under specified conditions, many drug substances can be identified by their R_f values and by the colors that their TLC spots turn when treated with various reagents or when observed under certain visualization methods. In the experiment that follows, TLC is applied to the analysis of an unknown analgesic drug.

REFERENCES

- Burtis, C. A.; and Bruns, D. E. *Tietz Fundamentals of Clinical Chemistry*, 4th ed.; W. B. Saunders: Philadelphia, 2014.
- James, S. H.; Nordby, J. J.; and Bell, S. *Forensic Science: An Introduction to Scientific and Investigative Techniques*, 4th ed.; CRC Press: Boca Raton, FL, 2014.
- Keller, E. Origin of Modern Criminology. *Chemistry* **1969**, *42*, 8.
- Keller, E. Forensic Toxicology: Poison Detection and Homicide. *Chemistry* **1970**, *43*, 14.
- Lieu, V. T. Analysis of APC Tablets. *J. Chem. Educ.* **1971**, *48*, 478.
- Neman, R. L. Thin Layer Chromatography of Drugs. *J. Chem. Educ.* **1972**, *49*, 834.
- Rodgers, S. S. Some Analytical Methods Used in Crime Laboratories. *Chemistry* **1969**, *42*, 29.
- A collection of articles on forensic chemistry can be found in
Berry, K.; and Outlaw, H. E., Eds. Forensic Chemistry—A Symposium Collection. *J. Chem. Educ.* **1985**, *62* (Dec), 1043–1065.

TLC Analysis of Analgesic Drugs

Thin-layer chromatography

In this experiment, thin-layer chromatography (TLC) will be used to determine the composition of various over-the-counter analgesics. If the instructor chooses, you may also be required to identify the components and actual identity (trade name) of an unknown analgesic. You will be given two commercially prepared TLC plates with a flexible backing and a silica-gel coating with a fluorescent indicator. On the first TLC plate, a reference plate, you will spot five standard compounds often used in analgesic formulations. In addition, a standard reference mixture containing four of these same compounds will be spotted. Ibuprofen is omitted from this standard mixture because it would overlap with salicylamide after the plate is developed. On the second plate (the sample plate), you will spot Naproxen sodium as an additional standard and four commercial analgesic preparations in order to determine their composition. At your instructor's option, one or more of these may be an unknown.

The standard compounds will all be available as solutions of 1 g of each dissolved in 20 mL of a 50:50 mixture of methylene chloride and ethanol. The purpose of the first reference plate is to determine the order of elution (R_f values) of the known substances and to index the standard reference mixture. Several of the substances have similar R_f values, but you will note a different behavior for each spot with the visualization methods. On the sample plate, the standard reference mixture will be spotted, along with Naproxen sodium and several solutions that you will prepare from commercial analgesic tablets. These tablets will each be crushed and dissolved in a 50:50 methylene chloride–ethanol mixture for spotting.

Reference Plate		Sample Plate	
Acetaminophen	(Ac)	Naproxen sodium	(Nap)
Aspirin	(Asp)	Sample 1*	(1)
Caffeine	(Cf)	Sample 2*	(2)
Ibuprofen	(Ibu)	Sample 3*	(3)
Salicylamide	(Sal)	Sample 4*	(4)
Reference mixture	(Ref)	Reference mixture	(Ref)

*At the instructor's option, one or more of the samples may be an unknown.

Two methods of visualization will be used to observe the positions of the spots on the developed TLC plates. First, the plates will be observed while under illumination from a short-wavelength ultraviolet (UV) lamp. This is best done in a darkened room or in a fume hood that has been darkened by taping butcher paper or aluminum foil over the lowered glass cover. Under these conditions, some of the spots will appear as dark areas on the plate, while others will fluoresce brightly. This difference in appearance under UV illumination will help to distinguish the

substances from one another. You will find it convenient to outline very lightly in *pencil* the spots observed and to place a small **x** inside those spots that fluoresce. For a second means of visualization, iodine vapor will be used. Not all the spots will become visible when treated with iodine, but some will develop yellow, tan, or deep brown colors. The differences in the behaviors of the various spots with iodine can be used to further differentiate among them.

It is possible to use several developing solvents for this experiment, but ethyl acetate with 0.5% glacial acetic acid added is preferred. The small amount of glacial acetic acid supplies protons and suppresses ionization of aspirin, ibuprofen, and naproxen sodium, allowing them to travel upward on the plates in their protonated form. Without the acid, these compounds do not move.

In some analgesics, you may find ingredients besides the five mentioned previously. Some include an antihistamine and some contain a mild sedative. For instance, Midol contains N-cinnamylephedrine (cinnamedrine), an antihistamine, and Excedrin PM contains the sedative methapyrilene hydrochloride. Cope contains the related sedative methapyrilene fumarate. Some tablets may be colored with a chemical dye.

REQUIRED READING

<i>Review:</i>	Essay	Analgesics
<i>New:</i>	Technique 19	Column Chromatography, Sections 19.1–19.3
	Technique 20	Thin-Layer Chromatography
	Essay	Identification of Drugs

SPECIAL INSTRUCTIONS

You must first examine the developed plates under ultraviolet light. After comparisons of *all* plates have been made with UV light, iodine vapor can be used. The iodine permanently affects some of the spots, making it impossible to go back and repeat the UV visualization. Take special care to notice those substances that have similar R_f values; these spots each have a different appearance when viewed under UV illumination or a different staining color with iodine, allowing you to distinguish among them.

Aspirin presents some special problems because it is present in a large amount in many of the analgesics and because it hydrolyzes easily. For these reasons, the aspirin spots often show excessive tailing.

SUGGESTED WASTE DISPOSAL

Dispose of all development solvent in the container for nonhalogenated organic solvents. Dispose of the ethanol–methylene chloride mixture in the container for halogenated organic solvents. The micropipettes used for spotting the solution should be placed in a special container labeled for that purpose. The TLC plates should be stapled in your lab notebook.

NOTES TO THE INSTRUCTOR

If you wish, students may work in pairs on this experiment, each one preparing one of the two plates.


Perform the thin-layer chromatography with flexible Silica Gel 60 F-254 plates (EM Science, No. 5554-7). If the TLC plates have not been purchased recently, you should place them in an oven at 100°C for 30 minutes and store them in a desiccator until used. If you use different thin-layer plates, try out the experiment before using them with a class. Other plates may not resolve all five substances.

Ibuprofen and salicylamide have approximately the same R_f value, but they show up differently under the detection methods. For reasons that are not yet clear, ibuprofen sometimes gives two or even three spots. Naproxen sodium has approximately the same R_f as aspirin. Once again, however, these analgesics show up differently under the detection methods. Fortunately, naproxen sodium is not combined with aspirin or ibuprofen in any current commercial product.

PROCEDURE

Initial Preparations

You will need at least 12 capillary micropipettes to spot the plates. The preparation of these pipettes is described and illustrated in Technique 20, Section 20.4. A common error is to pull the center section out too far when making these pipettes, with the result that too little sample is applied to the plate. If this happens, you won't see *any* spots. Follow the directions carefully.

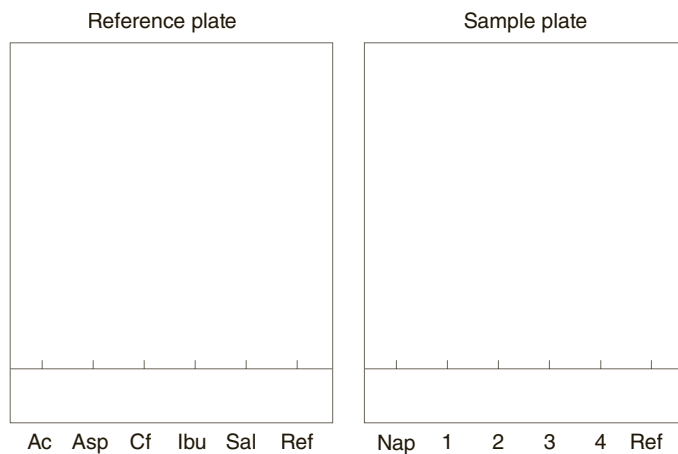

Preparing TLC Plates.

After preparing the micropipettes, obtain two 100-cm × 6.6-cm TLC plates (EM Science Silica Gel 60 F-254, No. 5554-7) from your instructor. These plates have a flexible backing, but they should not be bent excessively. Handle them carefully or the adsorbent may flake off. Also, you should handle them only by the edges; the surface should not be touched. Using a lead pencil (not a pen), *lightly* draw a line across the plates (short dimension) about 1 cm from the bottom. Using a centimeter ruler, move its index about 0.6 cm in from the edge of the plate and lightly mark off six 1-cm intervals on the line (see figure above). These are the points at which the samples will be spotted. If you are preparing two reference plates, it would be a good idea to mark a small number 1 or 2 in the upper right-hand corner of each plate to allow easy identification.

Spotting the First Reference Plate

On the first plate, starting from left to right, spot acetaminophen, then aspirin, caffeine, ibuprofen, and salicylamide. This order is alphabetic and will avoid any further memory problems or confusion. Solutions of these compounds will be found in

small bottles on the supply shelf. The standard reference mixture (Ref) also found on the supply shelf, is spotted in the last position. The correct method of spotting a TLC plate is described in Technique 20, Section 20.4. It is important that the spots be made as small as possible (ca. 1–2 mm in diameter). With too much sample, the spots will tail and will overlap one another after development. With too little sample, no spots will be observed after development. The optimum applied spot should be about 1–2 mm (1/6 in.) in diameter. If scrap pieces of the TLC plates are available, it would be a good idea to practice spotting on these before preparing the actual sample plates.

Preparing the Development Chamber

When the reference plate has been spotted, obtain a 16-oz wide-mouth, screw-cap jar (or other suitable container) for use as a development chamber. The preparation of a development chamber is described in Technique 20, Section 20.5. Because the backing on the TLC plates is very thin, if they touch the filter paper liner of the development chamber *at any point*, solvent will begin to diffuse onto the absorbent surface at that point. To avoid this, you may either omit the liner or make the following modification.

If you wish to use a liner, use a very narrow strip of filter paper (approximately 5 cm wide). Fold it into an L shape that is long enough to traverse the bottom of the jar and extend up the side to the top of the jar. TLC plates placed in the jar for development should *straddle* this liner strip, but not touch it.

When the development chamber has been prepared, obtain a small amount of the development solvent (0.5% glacial acetic acid in ethyl acetate). Your instructor should prepare this mixture; it contains such a small amount of acetic acid that small individual portions are difficult to prepare. Fill the chamber with the development solvent to a depth of about 0.5–0.7 cm. If you are using a liner, be sure it is saturated with the solvent. Recall that the solvent level must not be above the spots on the plate or the samples will dissolve off the plate into the reservoir instead of developing.

Development of the Reference TLC Plate

Place the spotted plate (or plates) in the chamber (straddling the liner if one is present) and allow the spots to develop. If you are doing two reference plates, both plates may be placed in the same development jar. Be sure the plates are placed in the developing jar so that their bottom edge is parallel to the bottom of the jar (straight, not tilted); if not, the solvent front will not advance evenly, increasing the difficulty of making good comparisons. The plates should face each other and slant or lean back in opposite directions. When the solvent has risen to a level about 0.5 cm from the top of the plate, remove each plate from the chamber (in the hood) and, using a lead pencil, mark the position of the solvent front. Set the plate on a piece of paper towel to dry. It may be helpful to place a small object under one end to allow optimum air flow around the drying plate.

UV Visualization of the Reference Plate

When the plate is dry, observe it under a short-wavelength UV lamp, preferably in a darkened hood or a darkened room. Lightly outline all of the observed spots with a pencil. Carefully notice any differences in behavior between the spotted substances. Several compounds have similar R_f values, but the spots have a different appearance under UV illumination or iodine staining. Currently, there are no commercial analgesic preparations containing any compounds that have the same R_f values, but you will need to be able to distinguish them from one another to identify which

one is present. Before proceeding, make a sketch of the plates in your notebook and note the differences in appearance that you observed. Using a ruler marked in millimeters, measure the distance that each spot has traveled relative to the solvent front. Calculate R_f values for each spot (see Technique 20, Section 20.9).

Analysis of Commercial Analgesics or Unknowns (Sample Plate)

Next, obtain half a tablet of each of the analgesics to be analyzed on the final TLC plate. If you were issued an unknown, you may analyze four other analgesics of your choice; if not, you may analyze five. The experiment will be most interesting if you make your choices in a way that gives a wide spectrum of results. Try to pick at least one analgesic each containing aspirin, acetaminophen, ibuprofen, a newer analgesic, and, if available, salicylamide. If you have a favorite analgesic, you may wish to include it among your samples. Take each analgesic half-tablet, place it on a smooth piece of notebook paper, and crush it well with a spatula. Transfer each crushed half-tablet to a labeled test tube or a small Erlenmeyer flask. Using a graduated cylinder, mix 15 mL of absolute ethanol and 15 mL of methylene chloride. Mix the solution well. Add 5 mL of this solvent to each of the crushed half-tablets and then heat each of them *gently* for a few minutes on a steam bath or sand bath at about 100°C. Not all of the tablet material will dissolve, because the analgesics usually contain an insoluble binder. In addition, many of them contain inorganic buffering agents or coatings that are insoluble in this solvent mixture. After heating the samples, allow them to settle and then spot the clear liquid extracts (1–4) on the sample plate. Spot the standard solution of naproxen on the left-hand edge, and spot the standard reference solution (Ref) on the right-hand edge of the plate (see figure above). Develop the plate in 0.5% glacial acetic acid–ethyl acetate as before. Observe the plate under UV illumination and mark the visible spots as you did for the first plate. Sketch the plate in your notebook and record your conclusions about the contents of each tablet. This can be done by directly comparing your plate to the reference plate(s)—they can all be placed under the UV light at the same time. If you were issued an unknown, try to determine its identity (trade name).

Iodine Analysis

Do not perform this step until UV comparisons of all the plates are complete. When ready, place the plates in a jar containing a few iodine crystals, cap the jar, and warm it gently on a steam bath or warm hot plate until the spots begin to appear. Notice which spots become visible and note their relative colors. You can directly compare colors of the reference spots to those on the unknown plate(s). Remove the plates from the jar and record your observations in your notebook.

QUESTIONS


1. What happens if the spots are made too large when preparing a TLC plate for development?
2. What happens if the spots are made too small when preparing a TLC plate for development?
3. Why must the spots be above the level of the development solvent in the developing chamber?
4. What would happen if the spotting line and positions were marked on the plate with a ball-point pen?
5. Is it possible to distinguish two spots that have the same R_f value but represent different compounds? Give two different methods.
6. Name some advantages of using acetaminophen (Tylenol) instead of aspirin as an analgesic.

Caffeine

The origins of coffee and tea as beverages are so old that they are lost in legend. Coffee is said to have been discovered by an Abyssinian goatherd who noticed an unusual friskiness in his goats when they consumed a certain little plant with red berries. He decided to try the berries himself and discovered coffee. The Arabs soon cultivated the coffee plant, and one of the earliest descriptions of its use is found in an Arabian medical book circa AD 900. The great systematic botanist Linnaeus named the plant *Coffea arabica*.

One legend of the discovery of tea—from the Orient, as you might expect—attributes the discovery to Daruma, the founder of Zen. Legend has it that he inadvertently fell asleep one day during his customary meditations. To be assured that this indiscretion would not recur, he cut off both eyelids. Where they fell to the ground, a new plant took root that had the power to keep a person awake. Although some experts assert that the medical use of tea was reported as early as 2737 BC in the pharmacopeia of Shen Nung, an emperor of China, the first indisputable reference is from the Chinese dictionary of Kuo P'o, which appeared in AD 350. The nonmedical, or popular, use of tea appears to have spread slowly. Not until about AD 700 was tea widely cultivated in China. Tea is native to upper Indochina and upper India, so it must have been cultivated in these places before its introduction to China. Linnaeus named the tea shrub *Thea sinensis*; however, tea is more properly a relative of the camellia, and botanists have renamed the shrub *Camellia thea*.

The active ingredient that makes tea and coffee valuable to humans is caffeine. Caffeine is an **alkaloid**, a class of naturally occurring compounds containing nitrogen and having the properties of an organic amine base (alkaline, hence, *alkaloid*). Tea and coffee are not the only plant sources of caffeine. Others include kola nuts, maté leaves, guarana seeds, and, in small amount, cocoa beans. The pure alkaloid was first isolated from coffee in 1821 by the French chemist Pierre Jean Robiquet.

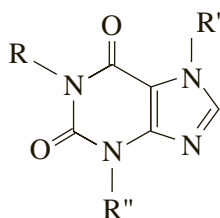

XANTHINES

Xanthine R = R' = R'' = H

Caffeine R = R' = R'' = CH₃

Theophylline R = R'' = CH₃, R' = H

Theobromine R = H, R' = R'' = CH₃

Caffeine belongs to a family of naturally occurring compounds called **xanthines**. The xanthines, in the form of their plant progenitors, are possibly the oldest known stimulants. They all, to varying extents, stimulate the central nervous system and the skeletal muscles. This stimulation results in an increased alertness, the ability to put off sleep, and an increased capacity for thinking. Caffeine is the most powerful xanthine in this respect. It is the main ingredient of the popular No-Doz keep-alert tablets. Although caffeine has a powerful effect on the central nervous system, not all xanthines are as effective. Thus, theobromine, the xanthine found in cocoa, has fewer central

nervous system effects. It is, however, a strong **diuretic** (induces urination) and is useful to doctors in treating patients with severe water-retention problems. Theophylline, a second xanthine found in tea, also has fewer central nervous system effects but is a strong **myocardial** (heart muscle) stimulant; it **dilates** (relaxes) the coronary artery that supplies blood to the heart. Its most important use is in the treatment of bronchial asthma, because it has the properties of a **bronchodilator** (relaxes the bronchioles of the lungs). Because it is also a **vasodilator** (relaxes blood vessels), it is often used in treating hypertensive headaches. It is also used to alleviate and to reduce the frequency of attacks of **angina pectoris** (severe chest pain). In addition, it is a more powerful diuretic than theobromine.

One can develop both a tolerance for the xanthines and a dependence on them, particularly caffeine. The dependence is real, and a heavy user (>5 cups of coffee per day) will experience lethargy, headache, and perhaps nausea after about 18 hours of abstinence. An excessive intake of caffeine may lead to restlessness, irritability, insomnia, and muscular tremor. Caffeine can be toxic, but to achieve a lethal dose of caffeine, one would have to drink about 100 cups of coffee over a relatively short period.

Caffeine is a natural constituent of coffee, tea, and kola nuts (*Kola nitida*). Theophylline is found as a minor constituent of tea. The chief constituent of cocoa is theobromine. The amount of caffeine in tea varies from 2% to 5%. In one analysis of black tea, the following compounds were found: caffeine, 2.5%; theobromine, 0.17%; theophylline, 0.013%; adenine, 0.014%; and guanine and xanthine, traces. Coffee beans can contain up to 5% by weight of caffeine, and cocoa contains around 5% theobromine. Commercial cola is a beverage based on a kola nut extract. We cannot easily get kola nuts in this country, but we can get the ubiquitous commercial extract as a syrup. The syrup can be converted into “cola.” The syrup contains caffeine, tannins, pigments, and sugar. Phosphoric acid is added, and caramel is added to give the syrup a deep color. The final drink is prepared by adding water and carbon dioxide under pressure, to give the bubbly mixture. Before decaffeination, the Food and Drug Administration required a “cola” to contain some caffeine (about 0.2 mg per ounce). In 1990, when new nutrition labels were adopted, this requirement was dropped. The Food and Drug Administration again currently requires that a “cola” contain *some* caffeine, but limits this amount, to a maximum of 5 milligrams per ounce. To achieve a regulated level of caffeine, most manufacturers remove all caffeine from the kola extract and then re-add the correct amount to the syrup. The caffeine content of various beverages is listed in the accompanying table.

Given the recent popularity of gourmet coffee beans and espresso stands, it is interesting to consider the caffeine content of these specialty beverages. Gourmet coffee certainly has more flavor than the typical ground coffee you may find on any grocery store shelf, and the concentration of brewed gourmet coffee tends to be higher than ordinary drip-grind coffee. Brewed gourmet coffee probably contains something on the order of 20–25 mg of caffeine per ounce of liquid. Espresso coffee is a very concentrated, dark-brewed coffee. Although the darker roasted beans used for espresso actually contain less caffeine per gram than regularly roasted beans, the method of preparing espresso (extraction using pressurized steam) is more efficient, and a higher percentage of the total caffeine in the beans is extracted. The caffeine content per ounce of liquid, therefore, is substantially higher than in most brewed coffees. The serving size for espresso coffee, however, is much smaller than for ordinary coffee (about 1.5–2 oz per serving), so the total caffeine available in a serving of espresso turns out to be about the same as in a serving of ordinary coffee.

Amount of Caffeine (mg/oz) Found in Beverages

Brewed coffee	12–30	Tea	4–20
Instant coffee	8–20	Cocoa (but 20 mg/oz theobromine)	0.5–2
Espresso (1 serving = 1.5–2 oz)	50–70	Coca-Cola	3.75
Decaffeinated coffee	0.4–1.0		

Note: The average cup of coffee or tea contains about 5–7 ounces of liquid. The average bottle of cola contains about 12 ounces of liquid.

Because of the central nervous system effects from caffeine, many people prefer **decaffeinated** coffee. The caffeine is removed from coffee by extracting the whole beans with an organic solvent. Then the solvent is drained off, and the beans are steamed to remove any residual solvent. The beans are dried and roasted to bring out the flavor. Decaffeination reduces the caffeine content of coffee to the range of 0.03% to 1.2%. The extracted caffeine is used in various pharmaceutical products, such as APC tablets.

Among coffee lovers there is some controversy about the best method to remove the caffeine from coffee beans. **Direct contact** decaffeination uses an organic solvent (usually methylene chloride) to remove the caffeine from the beans. When the beans are subsequently roasted at 200°C, virtually all traces of the solvent are removed, because methylene chloride boils at 40°C. The advantage of direct contact decaffeination is that the method removes only the caffeine (and some waxes), but leaves the substances responsible for the flavor of the coffee intact in the bean. A disadvantage of this method is that all organic solvents are toxic to some extent.

Water process decaffeination is favored among many drinkers of decaffeinated coffee because it does not use organic solvents. In this method, hot water and steam are used to remove caffeine and other soluble substances from the coffee. The resulting solution is then passed through activated charcoal filters to remove the caffeine. Although this method does not use organic solvents, the disadvantage is that water is not a very selective decaffeinating agent. Many of the flavor oils in the coffee are removed at the same time, resulting in a coffee with a somewhat bland flavor.

A third method, the **carbon dioxide decaffeination process**, is being used with increasing frequency. The raw coffee beans are moistened with steam and water, and they are then placed into an extractor where they are treated with carbon dioxide gas under very high temperature and pressure. Under these conditions, the carbon dioxide gas is in a **supercritical** state, which means that it takes on the characteristics of both a liquid and a gas. The supercritical carbon dioxide acts as a selective solvent for caffeine, thus extracting it from the beans.

There are, however, benefits to ingesting caffeine. Small amounts of caffeine have been found to be helpful in controlling weight, alleviating pain, and reducing the symptoms of asthma and other breathing problems. Recently, studies on mice indicate that caffeine may help to reverse or slow the development of Alzheimer's disease in mice. Other studies on humans indicate that caffeine may reduce the likelihood of developing Parkinson's disease and reduce the risk of colon cancer.

Another problem, posed by the beverage tea, is that in some cases persons who consume high quantities of tea may show symptoms of Vitamin B1 (thiamine) deficiency. It is suggested that the tannins in the tea may complex with the thiamine, rendering it unavailable for use. An alternative suggestion is that caffeine may reduce the levels of the enzyme transketolase, which depends on the presence of thiamine for its activity. Lowered levels of transketolase would produce the same symptoms as lowered levels of thiamine.

REFERENCES

- Arendash, G. W.; Mori, T.; Cao, C.; Mamcarz, M.; Runfeldt, M.; Dickson, A.; Rezai-Zadeh, K.; Tan, J.; Citron, B. A.; Lin, X.; et al. Caffeine Reverses Cognitive Impairment and Decreases Brain Amyloid- β Levels in Aged Alzheimer's Disease Mice. *J. Alzheimers Dis.* **2009**, *17*, 661–680.
- Emboden, W. *The Stimulants. Narcotic Plants*, rev. ed.; Macmillan: New York, 1979.
- Hart, C.; Ksir, C.; and Ray, O. Caffeine. *Drugs, Society, and Human Behavior*, 15th ed.; C. V. Mosby: St. Louis, 2012.
- Ritchie, J. M. Central Nervous System Stimulants. II: The Xanthines. In *The Pharmacological Basis of Therapeutics*, 8th ed.; Goodman L. S., Gilman, A., Eds.; Macmillan: New York, 1990.
- Ross, G. W.; Abbott, R. D.; Petrovich, H.; Morens, D. M.; Grandinetti, A.; Tung, K-H.; Tanner, C. M.; Masaki, K. H.; Blanchette, P. L.; Curb J. D.; et al. Association of Coffee and Caffeine Intake with the Risk of Parkinson Disease. *J. Am. Med. Assoc.* **2000**, *283* (May 24), 2674–2679.
- Taylor, N. *Plant Drugs That Changed the World*; Dodd Mead: New York, 1965; pp. 54–56.
- Taylor, N. Three Habit-Forming Nondangerous Beverages. *Narcotics—Nature's Dangerous Gifts*; Dell: New York, 1970. (Paperbound revision of Flight from Reality.)

Isolation of Caffeine from Tea or Coffee

Isolation of a natural product

Extraction

Sublimation

Green chemistry



Solid phase extraction (optional)


Gas chromatography (optional)

Infrared spectroscopy (optional)

In Experiment 13A, caffeine is isolated from tea leaves. The chief problem with the isolation is that caffeine does not exist alone in tea leaves but is accompanied by other natural substances from which it must be separated. The main component of tea leaves is cellulose, which is the principal structural material of all plant cells. Cellulose is a polymer of glucose. Because cellulose is virtually insoluble in water, it presents no problems in the isolation procedure. Caffeine, on the other hand, is water soluble and is one of the main substances extracted into the solution called tea. Caffeine constitutes as much as 5% by weight of the leaf material in tea plants.

Tannins also dissolve in the hot water used to extract tea leaves. The term **tannin** does not refer to a single homogeneous compound or even to substances that have similar chemical structure. It refers to a class of compounds that have certain properties in common. Tannins are phenolic compounds having molecular weights between 500 and 3000. They are widely used to tan leather. They precipitate alkaloids and proteins from aqueous solutions. Tannins are usually divided into two classes: those that can be **hydrolyzed** (react with water) and those that cannot. Tannins of the first type that are found in tea generally yield glucose and gallic acid when they are hydrolyzed. These tannins are esters of gallic acid and glucose. They represent structures in which some of the hydroxyl groups in glucose have been esterified by digalloyl groups. The nonhydrolyzable tannins found in tea are condensation polymers of catechin. These polymers are not uniform in structure; catechin molecules are usually linked at ring positions 4 and 8.

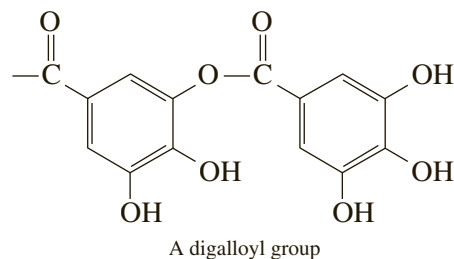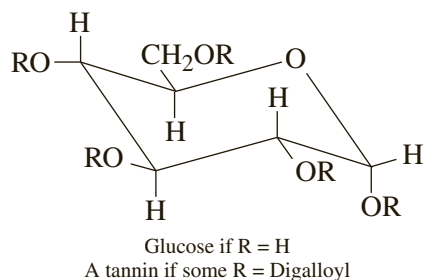


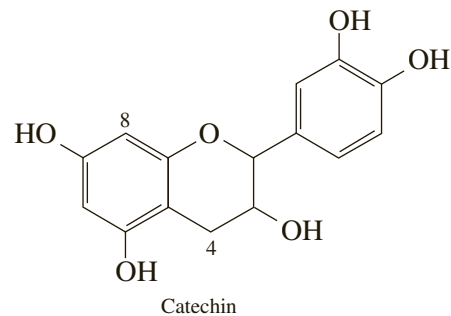

When tannins are extracted into hot water, some of these compounds are partially hydrolyzed to form free gallic acid. The tannins, because of their phenolic groups, and gallic acid, because of its carboxyl groups, are both acidic. If sodium carbonate, a base, is added to tea water, these acids are converted to their sodium salts that are highly soluble in water.

Although caffeine is soluble in water, it is much more soluble in the organic solvent methylene chloride. Caffeine can be extracted from the basic tea solution with methylene chloride, but the sodium salts of gallic acid and the tannins remain in the aqueous layer.

The brown color of a tea solution is due to flavonoid pigments and chlorophylls and to their respective oxidation products. Although chlorophylls are soluble in methylene chloride, most other substances in tea are not. Thus, the methylene chloride extraction of the basic tea solution removes nearly pure caffeine. The methylene chloride is easily removed by evaporation (bp 40°C) to leave the crude caffeine. The caffeine is then purified by sublimation at reduced pressure to prevent decomposition.

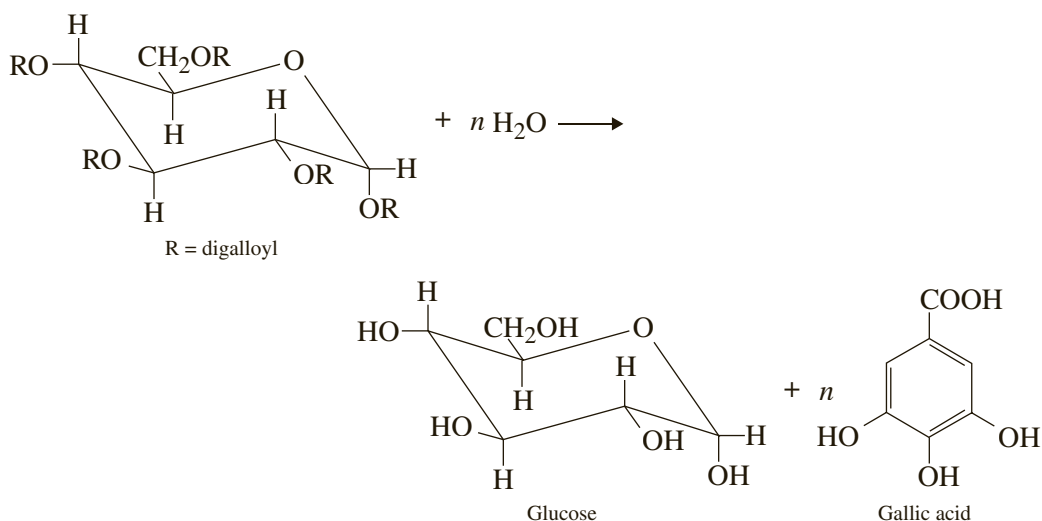

The methylene chloride solvent used in Experiment 13A is toxic and a suspected carcinogen, but it remains as one of the best solvents for extracting caffeine from aqueous tea and coffee solutions. Other solvents simply do not do as good a job of extracting caffeine. Besides being toxic, chlorinated solvents are also environmental pollutants (see the “Green Chemistry” essay). This solvent must be disposed of correctly to avoid environmental problems. Ethyl acetate is “more

Green" than methylene chloride, but unfortunately it does not extract caffeine from aqueous solutions nearly as efficiently as methylene chloride.

Solid phase extraction (SPE), described in Technique 12, Section 12.14, avoids the use of methylene chloride, and the technique will be employed in Experiment 13B for extracting caffeine from tea or coffee. The material used in this SPE procedure is prepared by alkylating silica with 18-carbon chains. This process changes the silica from a polar substance to a relatively nonpolar one. The resulting material, called C-18 silica, is placed in a tube that resembles a syringe body. Because the silica is now nonpolar, water and other polar substances such as tannins and gallic acid will pass through the column, whereas the C-18 silica retains caffeine. To remove the caffeine from the column, ethyl acetate is added.

REQUIRED READING

<i>Review:</i>	Techniques 5 and 6	
	Technique 7	Reaction Methods, Section 7.10
	Technique 9	Physical Constants of Solids: The Melting Point
<i>New:</i>	Technique 12	Extractions, Separations, and Drying Agents
	Technique 17	Sublimation
	Essay	Caffeine
	Essay	Green Chemistry
	Technique 22	Gas Chromatography (optional)
	Technique 25	Infrared Spectroscopy (optional)
	Technique 28	Mass Spectrometry (optional)

SPECIAL INSTRUCTIONS

Be careful when handling methylene chloride. It is a toxic solvent and suspected carcinogen. You should not breathe it or spill it on yourself. Wear gloves and work in a fume hood. It is recommended that you use glass centrifuge tubes with screw caps. Check them first with water to make sure that they will not leak while you shake them. Plastic centrifuge tubes with screw caps can be used instead of glass ones. However, there is some tendency for the plastic to dissolve in methylene chloride, and they should not be in contact with methylene chloride for long periods (more than 1 hour).

SUGGESTED WASTE DISPOSAL

You must dispose of methylene chloride in a waste container marked for the disposal of halogenated organic waste. Dispose of the tea bag in a trash can. The aqueous solutions obtained after the extraction steps must be disposed of in a waste container labeled for aqueous waste.

13A EXPERIMENT 13A

Extraction of Caffeine from Tea with Methylene Chloride

PROCEDURE

Preparing the Tea Solution

Place 20 mL of water in a 50-mL beaker. Cover the beaker with a watch glass and heat the water on a hot plate until it is almost boiling. From information provided by your instructor, determine the weight of the tea present in the tea bag.¹ Place the tea bag into the hot water so that it lies flat on the bottom of the beaker and is covered as completely as possible with water. Replace the watch glass and continue heating for about 15 minutes. During this heating period, it is important to push down *gently* on the tea bag with a test tube so that all the tea leaves are in constant contact with water. As the water evaporates during this heating step, replace it by adding water from a Pasteur pipette.

Using the Pasteur pipette, transfer the concentrated tea solution to two centrifuge tubes fitted with screw caps. Try to keep the liquid volume in each centrifuge tube approximately equal. To squeeze additional liquid out of the tea bag, hold the tea bag on the inside wall of the beaker and roll a test tube back and forth while exerting *gentle* pressure on the tea bag. Press out as much liquid as possible without breaking the bag. Combine this liquid with the solution in the centrifuge tubes. Place the tea bag on the bottom of the beaker again and pour 2 mL of hot water over the bag. Squeeze the liquid out, as just described, and transfer this liquid to the centrifuge tubes. Add 0.5 g of sodium carbonate to the hot liquid in each centrifuge tube. Cap the tubes and shake the mixture until the solid dissolves.

Extraction and Drying

Cool the tea solution to room temperature. Using a calibrated Pasteur pipette, add 3 mL of methylene chloride to each centrifuge tube to extract the caffeine (Technique 12, Section 12.4). Cap the centrifuge tubes and gently shake the mixture for several seconds. Vent the tubes to release the pressure, being careful that the liquid does not squirt out toward you. Shake the mixture for an additional 30 seconds with occasional venting. To separate the layers and break the emulsion (see Technique 12, Section 12.10), centrifuge the mixture for several minutes (be sure to balance the centrifuge by placing the two centrifuge tubes on opposite sides). If an emulsion still remains (indicated by a green-brown layer between the clear methylene chloride layer and the top aqueous layer), centrifuge the mixture again.

Remove the lower organic layer with a Pasteur pipette and transfer it to a dry 25-mL Erlenmeyer flask. Be sure to squeeze the bulb before placing the tip of the Pasteur pipette into the liquid and try not to transfer any of the dark aqueous

¹ The weight of the tea in the bag varies, whereas the weight of the bag, string, and tag are relatively constant. Your instructor may provide you with the weight of everything except the tea itself. You can then determine the weight of the tea by directly weighing the tea bag provided you and subtracting the weight of the empty bag, string, and tag to determine the weight of the tea in the bag.

solution along with the methylene chloride layer. Add a fresh 3-mL portion of methylene chloride to the aqueous layer remaining in each centrifuge tube, cap the centrifuge tubes, and shake the mixture in order to carry out a second extraction. Separate the layers by centrifugation, as described previously. Combine the organic layers from each extraction into the flask containing the first extraction. If there are visible drops of the dark aqueous solution in the flask, transfer the methylene chloride solution to another flask using a clean, dry Pasteur pipette. If necessary, leave a small amount of the methylene chloride solution behind in order to avoid transferring any of the aqueous mixture. Add granular anhydrous sodium sulfate to dry the organic layer (Technique 12, Section 12.9). If all the sodium sulfate clumps together when the mixture is stirred with a spatula, add some additional drying agent. Allow the mixture to stand for 10–15 minutes. Stir occasionally with a spatula.

Evaporation

Transfer the dry methylene chloride solution with a Pasteur pipette to a dry, pre-weighed 25-mL Erlenmeyer flask, while leaving the drying agent behind. Evaporate the methylene chloride by heating the flask in a hot water bath (Technique 7, Section 7.10). This should be done in a hood and can be accomplished more rapidly if a stream of dry air or nitrogen gas is directed at the surface of the liquid. When the solvent is evaporated, the crude caffeine will coat the bottom of the flask. Do not heat the flask after the solvent has evaporated or you may sublime some of the caffeine. Weigh the flask and determine the weight of crude caffeine. Calculate the weight percentage recovery (see Technique 2, Section 2.2C) of caffeine from tea leaves, using the weight of tea that you started with. You may store the caffeine by simply placing a stopper firmly into the flask.

Sublimation of Caffeine

Caffeine can be purified by sublimation (Technique 17, Section 17.5). Assemble a sublimation apparatus as shown in Figure 17.2A.² Add approximately 0.5 mL of methylene chloride to the Erlenmeyer flask and transfer the solution to a clean, 5-mL, thin-walled, conical vial, using a clean and dry Pasteur pipette. Add a few more drops of methylene chloride to the flask in order to rinse the caffeine out completely. Transfer this liquid to the conical vial. Evaporate the methylene chloride from the conical vial by gentle heating in a warm water bath under a stream of dry air or nitrogen.

Insert the cold finger into the sublimation apparatus. If you are using the sublimator with the multipurpose adapter, adjust it so that the tip of the cold finger will be positioned about 1 cm above the bottom of the conical vial. Be sure that the inside of the assembled apparatus is clean and dry. If you are using an aspirator, install a trap between the aspirator and the sublimation apparatus. Turn on the vacuum and check to make sure that all joints in the apparatus are sealed tightly. Place *ice-cold* water in the inner tube of the apparatus. Heat the sample gently and carefully with a microburner to sublime the caffeine. Hold the burner in your hand (hold it at its base, *not* by the hot barrel), and apply the heat by moving the flame back and forth under the conical vial and up the sides. If the sample begins to melt, remove the flame for a few seconds before you resume heating. When sublimation is complete, discontinue heating. Remove the cold water and remaining ice from the inner tube and allow the apparatus to cool while continuing to apply the vacuum.

² If you are using another type of sublimation apparatus, your instructor will provide you with specific instructions on how to assemble it correctly.

When the apparatus is at room temperature, remove the vacuum and *carefully* remove the inner tube. If this operation is done carelessly, the sublimed crystals may be dislodged from the inner tube and fall back into the conical vial. Scrape the sublimed caffeine onto a tared piece of smooth paper and determine the weight of caffeine recovered. Calculate the weight percentage recovery (see Technique 2, Section 2.2C) of caffeine after the sublimation. Compare this value to the percentage recovery determined after the evaporation step. Determine the melting point of the purified caffeine. The melting point of pure caffeine is 236°C; however, the observed melting point will be lower. Submit the sample to the instructor in a labeled vial.

13B EXPERIMENT 13B (OPTIONAL)

Extraction of Caffeine from Tea or Coffee Using Solid Phase Extraction (SPE)

In this experiment, caffeine is isolated from brewed tea or coffee using solid phase extraction (Technique 12, Section 12.14). It is offered as a Green Chemistry alternative to the procedure given in Experiment 13A in which methylene chloride is used for the extraction of caffeine from tea.

The solid phase extraction (SPE) columns may be obtained commercially.³ The manufacturers covalently bond nonpolar alkyl groups to silica. This process converts the polar silica to a relatively nonpolar material. Often, 18-carbon alkyl groups are bonded to the silica creating a material that is referred to as C-18 silica.

In “normal” chromatography using nontreated silica, one expects polar substances to be attracted to the polar surface and to move more slowly through the column compared to nonpolar substances. This more traditional form of chromatography is referred to as **normal phase chromatography**. With C-18-treated silica, however, the polar substances will come off first, and the relatively nonpolar substance will be retained. This type of chromatography is referred to as **reverse phase chromatography**.

Experiment 12B makes use of reverse phase chromatography to remove caffeine from a tea or coffee solution. The technique is simple. See Technique 12, Figure 12.14, for a diagram of the apparatus that you will use. The procedure involves the following steps:

1. Condition the SPE reverse phase column with methanol and water.
2. Use a vacuum to draw an aqueous solution of brewed tea or coffee through the SPE reverse phase column.
3. Caffeine is retained by the solid phase column.
4. Water, tannins, gallic acid, and other polar substances pass through the column.
5. Remove the caffeine from the column with ethyl acetate.
6. Evaporate the ethyl acetate to yield the crude caffeine.
7. Purify the caffeine by sublimation.

³ We recommend a 6-mL column with 1000 mg of C18E sorbent (No. 8B-S001-JCH-S) for Experiment 13B. These columns are available from Phenomenex, 411 Madrid Avenue, Torrance, CA 90501-1430, phone: (310) 212-0555. About 25 mg of caffeine can be isolated from one batch of tea or coffee solution.

SUGGESTED WASTE DISPOSAL

This experiment approaches closely a high Green standard. Get into the spirit of the Green movement! Throw the tea bag into a trash can. Pour the aqueous solution from the column down the drain in a sink. Any remaining ethyl acetate should be given to another student who needs it or, as a last resort, poured into a nonhalogenated waste container. Remove the traces of ethyl acetate from the sodium sulfate drying agent with a blast of air and then dissolve the remaining solid inorganic material with water and pour it into an aqueous waste container (or pour it down the sink drain). Place the SPE column into a beaker that has been supplied for your use. Eventually, the tubes will be emptied and reused.

PROCEDURE

Preparing the Tea or Coffee Solution

Prepare a tea solution as described in Experiment 13A (first paragraph) and footnote 1. (If you desire to extract caffeine from coffee, dissolve about 1.5 g of instant coffee in 20 mL of boiling water.⁴) Using the Pasteur pipette, transfer the concentrated tea solution into two centrifuge tubes fitted with screw caps. Try to keep the liquid volume in each centrifuge tube approximately equal. To squeeze additional liquid out of the tea bag, hold the tea bag on the inside wall of the beaker and roll a test tube back and forth while exerting gentle pressure on the tea bag. Press out as much liquid as possible without breaking the bag. Combine this liquid with the solution in the centrifuge tubes. Place the tea bag on the bottom of the beaker again and pour 2 mL of hot water over the tea bag. Squeeze the liquid out as just described and transfer this liquid to the centrifuge tubes. Let the solution in the centrifuge tube cool to room temperature. When the solution is room temperature, cap the centrifuge tubes and centrifuge the samples for 4 to 5 minutes. When the centrifugation step is complete, a solid will be present at the bottom of the tubes. Transfer the solution with a Pasteur pipette into a small beaker so that the thick solid is left behind in the centrifuge tubes.

Preparing the SPE Column

The apparatus is shown in Technique 12, Figure 12.14. Insert a No. 1 neoprene adapter inside a No. 2 neoprene adapter and place both adapters in the top of a 50-mL filter flask. Clamp the flask securely and then apply a vacuum source. Place the SPE column⁵ inside the No. 1 neoprene adapter and condition the column under vacuum in the following way. Add 2 mL of methanol to the column, 1 mL at a time. Wait for the methanol to draw through the column completely before adding the second mL. Next, add 2 mL of room-temperature de-ionized (DI) water, 1 mL at a time. Proceed to the next step as quickly as possible. Don't let the column dry out.

Filtering the Tea Solution

While still drawing a vacuum, transfer the tea (or coffee) solution in the small beaker to the SPE column in approximately 1-mL portions, using a Pasteur pipette. If some solid is present in the beaker, be careful to avoid drawing up this solid at the bottom of the beaker (this solid will plug the SPE column). The SPE column should

⁴ Tea yields better-quality caffeine than that obtained from coffee. Caffeine from tea is relatively colorless, whereas the caffeine extracted from coffee is highly colored. About 25 mg of caffeine is isolated in either case. Sublimation removes much of the color from the tea and coffee samples.

⁵ See footnote 3 for the SPE column required for this experiment.

still be under vacuum. When all the solution has been drawn through the column, add 3 mL of room-temperature DI water in 1-mL portions. Continue applying the vacuum for 1 minute to remove as much water from the column as possible. Much of the caffeine is retained in the SPE column.

Removing the Caffeine from the SPE Column

Turn off the vacuum and carefully and slowly remove the hose from the side arm of the filter flask so that the liquid in the flask will not bump up onto the adapters and SPE column. Discard the filtrate in the filter flask. Clean and dry the filter flask that you used. Reconnect the SPE column to the apparatus and reapply the vacuum. Add a total of 7 mL of ethyl acetate to the column in approximately 1-mL portions. When all the ethyl acetate has been added, turn off the vacuum source and carefully remove the hose from the filter flask. Place the SPE column in the beaker supplied by the instructor.⁶ Add 1.6 g of granular anhydrous sodium sulfate to the ethyl acetate solution in the filter flask. This drying operation removes any remaining water from the ethyl acetate solution. Allow the mixture to stand for 15 minutes. Stir the solution occasionally with a spatula to complete the drying step (Technique 12, Section 12.9).

Isolation of the Caffeine

Transfer the dry ethyl acetate solution with a Pasteur pipette to a *dry, preweighed* 25-mL Erlenmeyer flask while leaving the drying agent behind. Rinse the drying agent with about 1 mL of fresh ethyl acetate and transfer it to the flask with your Pasteur pipette. Evaporate the ethyl acetate by heating the solution on a hot plate. This must be done in a hood and can be accomplished more rapidly if a stream of dry air or nitrogen gas is *carefully* directed at the surface of the liquid. Take the flask off the hot plate while some liquid remains and then use a gentle stream of air to remove the remaining traces of ethyl acetate from your sample. When the solvent has evaporated, crude solid caffeine will coat the bottom of the flask. Reweigh the flask to determine the weight of crude caffeine. Calculate the weight percentage recovery of caffeine obtained from your tea or coffee sample.⁷ You may store the caffeine by placing a stopper firmly into the flask.

Sublimation of Caffeine

The caffeine obtained from tea or coffee can be purified by sublimation using the procedure described in Experiment 13A. At your instructor's option, you may combine your sample with another student's sample for sublimation. After sublimation, determine the weight of caffeine recovered and calculate the weight percentage recovery of the caffeine. Compare this value to the amount of crude sample obtained. At the instructor's option, determine the melting point of the purified caffeine. The melting point of pure caffeine is 236°C; however, the observed melting point will be lower. Submit the sample to the instructor in a labeled vial unless it is to be used for infrared spectroscopy (recommended) or mass spectroscopy (also recommended).

Analysis of Caffeine

At the option of your instructor, obtain the infrared spectrum of the caffeine as a dry film by dissolving a sample in *dry* acetone (Technique 25, Section 25.4). Alternatively, a KBr pellet of the caffeine can be prepared to obtain the spectrum (Technique 25, Section 25.5). Include the infrared spectrum with your laboratory report, along with an interpretation of the principal peaks. A spectrum is included for comparison purposes.


⁶ *Note to the instructor:* It is recommended that the C-18 silica in the SPE columns not be reused.

⁷ See footnote 3 for information on the expected yield of caffeine from tea and coffee.

At the instructor's option, you may determine the mass spectrum of the caffeine by gas chromatography/mass spectrometry (GC-MS) (Technique 28). At the same time, you can determine the purity of your sample.⁸

Report

Attach your infrared spectra to your report and label the major peaks. If you determined the mass spectrum, try to identify the important fragment ion peaks (Technique 28). Include the melting point, if it was required. Report the weight percentage of the caffeine recovered from the tea or coffee sample before and after sublimation.

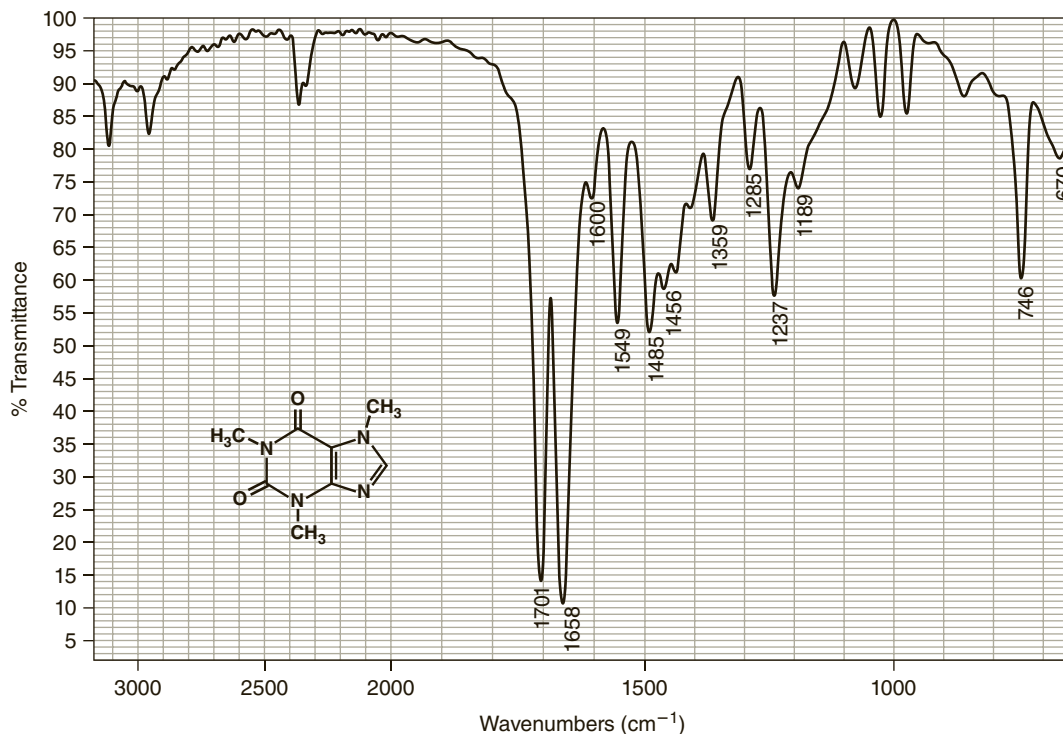

Spectrum of caffeine in KBr.


QUESTIONS

- Outline a separation scheme for isolating caffeine from tea. Use a flowchart similar in format to that shown in Technique 2, Section 2.2.
- Why was the sodium carbonate added in Experiment 13A?
- The crude caffeine isolated from tea often has a green tinge. Why?
- What are some possible explanations for why the melting point of your isolated caffeine was lower than the literature value (236°C)?
- An alternative procedure for removing the tannins and gallic acid is to heat the tea leaves in an aqueous mixture containing calcium carbonate. Calcium carbonate reacts with the tannins and gallic acid to form insoluble calcium salts of these acids. If this procedure were used, what additional step (not done in this experiment) would be needed in order to obtain an aqueous tea solution?
- What would happen to the caffeine if the sublimation step were performed at atmospheric pressure?

⁸The GC/MS may be obtained on the caffeine sample using a Varian 3800 GC with Saturn 2000 MS equipped with a capillary column that is 30 m long and 0.25 mm ID containing a 0.25- μ m film of Varian CP-Sil 8CB. The injector was set at 275°C, with He gas flow at 1 mL/min. The temperature of the oven was increased from 75 to 280°C at 20%/min and held at 280°C for 1 minute. The retention time for the caffeine is about 10 minutes under these conditions. No impurities were detected.

Esters—Flavors and Fragrances

Esters are a class of compounds widely distributed in nature. They have the general formula

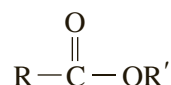











The simple esters tend to have pleasant odors. In many cases, although not exclusively so, the characteristic flavors and fragrances of flowers and fruits are due to compounds with the ester functional group. An exception is the case of the essential oils. The **organoleptic** qualities (odors and flavors) of fruits and flowers may often be due to a single ester, but more often, the flavor or the aroma is due to a complex mixture in which a single ester predominates. Some common flavor principles are listed in Table 1. Food and beverage manufacturers are familiar with these esters and often use them as additives to spruce up the flavor or odor of a dessert or beverage. Many times, such flavors or odors do not even have a natural basis, as is the case with the “juicy fruit” principle, isopentenyl acetate. An instant pudding that has the flavor of rum may never have seen its alcoholic namesake; this flavor can be duplicated by the proper admixture, along with other minor components, of ethyl formate and isobutyl propionate. The natural flavor and odor are not exactly duplicated, but most people can be fooled. Often, only a trained person with a high degree of gustatory perception, a professional taster, can tell the difference.

A single compound is rarely used in good-quality imitation flavoring agents. A formula for an imitation pineapple flavor that might fool an expert is listed in Table 2. The formula includes 10 esters and carboxylic acids that can easily be synthesized in the laboratory. The remaining seven oils are isolated from natural sources.

Flavor is a combination of taste, sensation, and odor transmitted by receptors in the mouth (taste buds) and nose (olfactory receptors). The stereochemical theory of odor is discussed in the essay that precedes Experiment 16. The four basic tastes (sweet, sour, salty, and bitter) are perceived in specific areas of the tongue. The sides of the tongue perceive sour and salty tastes, the tip is most sensitive to sweet tastes, and the back of the tongue detects bitter tastes. The perception of flavor, however, is not so simple. If it were, it would require only the formulation of various combinations of four basic substances—a bitter substance (a base), a sour substance (an acid), a salty substance (sodium chloride), and a sweet substance (sugar)—to duplicate any flavor! In fact, we cannot duplicate flavors in this way. The human possesses 9,000 taste buds. The combined response of these taste buds is what allows perception of a particular flavor.

Although the “fruity” tastes and odors of esters are pleasant, they are seldom used in perfumes or scents that are applied to the body. The reason for this is chemical. The ester group is not as stable under perspiration as the ingredients

Table 1 Ester Flavors and Fragrances

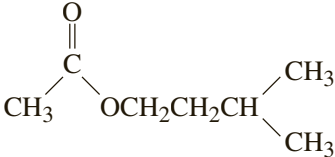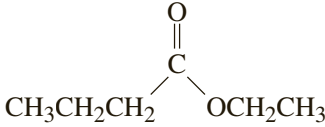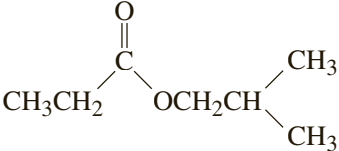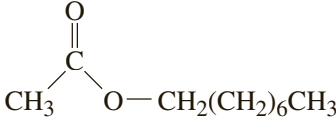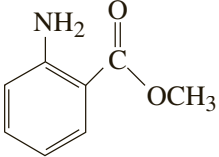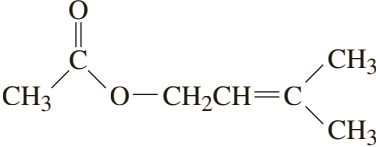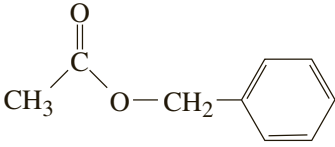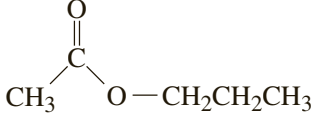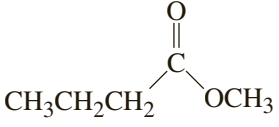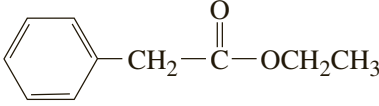 <p>Isoamyl acetate (banana) (alarm pheromone of honeybee)</p>	 <p>Ethyl butyrate (pineapple)</p>
 <p>Isobutyl propionate (rum)</p>	 <p>Octyl acetate (oranges)</p>
 <p>Methyl anthranilate (grape)</p>	 <p>Isopentenyl acetate ("Juicy Fruit")</p>
 <p>Benzyl acetate (peach)</p>	 <p>n-Propyl acetate (pear)</p>
 <p>Methyl butyrate (apple)</p>	 <p>Ethyl phenylacetate (honey)</p>

of the more expensive essential-oil perfumes. The latter are usually hydrocarbons (terpenes), ketones, and ethers extracted from natural sources. Esters, however, are used only for the cheapest toilet waters because on contact with sweat they undergo hydrolysis, giving organic acids. These acids, unlike their precursor esters, generally do not have a pleasant odor.

Butyric acid, for instance, has a strong odor like that of rancid butter (of which it is an ingredient) and is a component of what we normally call body odor. It is this substance that makes foul-smelling humans so easy for an animal to detect when downwind of them. It is also of great help to the bloodhound, which is trained to follow small traces of this odor.



Table 2 Artificial pineapple flavor

<i>Pure Compounds</i>	%	<i>Essential Oils</i>	%
Allyl caproate	5	Oil of sweet birch	1
Isoamyl acetate	3	Oil of spruce	2
Isoamyl isovalerate	3	Balsam Peru	4
Ethyl acetate	15	Volatile mustard oil	1
Ethyl butyrate	22	Oil cognac	5
Terpinyl propionate	3	Concentrated orange oil	4
Ethyl crotonate	5	Distilled oil of lime	<u>2</u>
Caproic acid	8		19
Butyric acid	12		
Acetic acid	<u>5</u>		
	81		

Ethyl butyrate and methyl butyrate, however, which are the *esters* of butyric acid, smell like pineapple and apple, respectively.

A sweet, fruity odor also has the disadvantage of possibly attracting fruit flies and other insects in search of food. Isoamyl acetate, the familiar solvent called banana oil, is particularly interesting. It is identical to a component of the alarm **pheromone** of the honeybee. Pheromone is the name applied to a chemical secreted by an organism that evokes a specific response in another member of the same species. This kind of communication is common among insects who otherwise lack means of intercourse. When a honeybee worker stings an intruder, an alarm pheromone, composed partly of isoamyl acetate, is secreted along with the sting venom. This chemical causes aggressive attack on the intruder by other bees, who swarm after the intruder. Obviously, it wouldn't be wise to wear a perfume compounded of isoamyl acetate near a beehive. Pheromones are discussed in more detail in the essay preceding Experiment 50.

REFERENCES

- The Givaudan Index*. Givaudan-Delawanna: New York, 1949. (Gives specifications of synthetics and isolates for perfumery.)
- Gould, R. F., Ed. *Flavor Chemistry, Advances in Chemistry*, Series No. 56; American Chemical Society: Washington, DC, 1966.
- Layman, P. L. Flavors and Fragrances Industry Taking on New Look. *Chem. Eng. News* **1987**, (Jul 20), 35.
- Moyler, D. Natural Ingredients for Flavours and Fragrances. *Chem. Ind.* **1991**, (Jan 7), 11.
- Rasmussen, P. W. Qualitative Analysis by Gas Chromatography—G.C. versus the Nose in Formulation of Artificial Fruit Flavors. *J. Chem. Edu.*, **1984**, 61 (Jan), 62.
- Shreve, R. N.; and Brink, J. *Chemical Process Industries*, 6th ed.; McGraw-Hill: New York, 1985. (A classic, now out of print, but probably still available in your library.)
- Surburg, H.; and Panten, J. *Common Fragrance and Flavor Materials*, 5th ed.; VCH Publishers: Weinheim, Germany, 2006.
- Welsh, F. W.; and Williams, R. E. Lipase Mediated Production of Flavor and Fragrance Esters from Fusel Oil. *J. Food Sci.*, **1989**, 54 (Nov/Dec), 1565.

Isopentyl Acetate (Banana Oil)

Esterification


Heating under reflux

Extraction

Simple distillation

Microscale boiling point

In this experiment you will prepare an ester, isopentyl acetate. This ester is often referred to as banana oil because it has the familiar odor of this fruit.

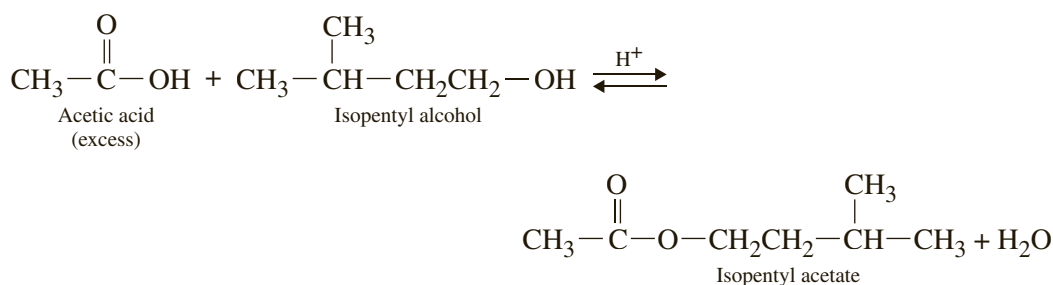

Isopentyl acetate is prepared by the direct esterification of acetic acid with isopentyl alcohol. Because the equilibrium does not favor the formation of the ester, it must be shifted to the right, in favor of the product, by using an excess of one of the starting materials. Acetic acid is used in excess because it is less expensive than isopentyl alcohol and more easily removed from the reaction mixture.

In the isolation procedure, much of the excess acetic acid and the remaining isopentyl alcohol are removed by extraction with sodium bicarbonate and water. After drying with anhydrous sodium sulfate, the ester is purified by distillation. The purity of the liquid product is analyzed by performing a microscale boiling point determination or infrared spectroscopy.

REQUIRED READING

<i>Review:</i>	Experiment 1	Introduction to Microscale Laboratory
	Techniques 5 and 6	
<i>New:</i>	Technique 7	Reaction Methods, Sections 7.2–7.4 and 7.6
	Technique 13	Physical Constants, Boiling Points
	Technique 12	Extractions, Separations, and Drying Agents
	Technique 14	Simple Distillation
	Essay	Esters—Flavors and Fragrances

If performing the optional infrared spectroscopy, also read:
Technique 25 Preparation of Samples for Spectroscopy

SPECIAL INSTRUCTIONS

Be careful when dispensing sulfuric and glacial acetic acids. They are corrosive and will attack your skin if you make contact with them. If you get one of these acids on your skin, wash the affected area with copious quantities of running water for 10–15 minutes.

Because a 1-hour reflux is required, you should start the experiment at the beginning of the laboratory period. During the reflux period, you may perform other work.

SUGGESTED WASTE DISPOSAL

Any aqueous solutions should be placed in a container specially designated for dilute aqueous waste. Place any excess ester in the nonhalogenated organic waste container.

NOTES TO THE INSTRUCTOR

Choose either Experiment 14A or Experiment 14B, but not both. The semimicroscale procedure requires the use of equipment not found in the typical microscale kit: a 20-mL round-bottom flask, a distillation head, and a vacuum takeoff adapter. The purpose of Experiment 14B is to allow an alternative to the use of a Hickman head for the distillation step.

This experiment has been carried out successfully using Dowex 50 × 2-100 ion-exchange resin instead of the sulfuric acid. Amberlyst-15 resin will also work.

14A EXPERIMENT 14A

Isopentyl Acetate (Microscale Procedure)

PROCEDURE

Apparatus

Using a 5-mL conical vial, assemble a reflux apparatus using a water-cooled condenser (Technique 7, Figure 7.2A). Top the condenser with a drying tube (Technique 7, Figure 7.10B) that contains a loose plug of glass wool. The purpose of the drying tube is to control odors rather than to protect the reaction from water. Use a hot plate and an aluminum block for heating.

Preparation

Remove the empty 5-mL conical vial, weigh it, and record its weight. Place approximately 1.0 mL of isopentyl alcohol ($MW = 88.2$, $d = 0.813$ g/mL) in the vial using an automatic pipette or a dispensing pump. Reweigh the vial containing the alcohol and subtract the tare weight to obtain an accurate weight for the alcohol. Add 1.5 mL of glacial acetic acid ($MW = 60.1$, $d = 1.06$ g/mL) using

an automatic pipette or dispensing pump. Using a disposable Pasteur pipette, add two to three drops of concentrated sulfuric acid. Swirl the liquid to mix. Add a small boiling stone (or a magnetic spin vane) and reattach the vial to the apparatus.

Reflux

Bring the mixture to a boil (aluminum block at about 150–160°C). Be sure to stir the mixture if you are using a spin vane instead of a boiling stone. Continue heating under reflux for 60–75 minutes. Remove the heating source and allow the mixture to cool to room temperature.

Workup

Disassemble the apparatus and, using a forceps, remove the boiling stone (or spin vane). Using a calibrated Pasteur pipette, slowly add 1.0 mL of 5% aqueous sodium bicarbonate to the cooled mixture in the conical vial. Stir the mixture in the vial with a microspatula until carbon dioxide evolution is no longer vigorous. Then cap the vial and shake *gently* with venting until the evolution of gas is complete. Using a Pasteur pipette, remove the lower aqueous layer and discard it. Repeat the extraction two more times, as outlined previously, using a fresh 1.0-mL portion of 5% sodium bicarbonate solution each time.

If droplets of water are evident in the vial containing the ester, transfer the ester to a dry conical vial using a dry Pasteur pipette. Dry the ester over granular anhydrous sodium sulfate (see Technique 12, Section 12.9). Allow the capped solution to stand for 10–15 minutes. Transfer the dry ester with a Pasteur pipette into a 3-mL conical vial while leaving the drying agent behind. If necessary, pick out any granules of sodium sulfate with the end of a spatula.

Distillation

Add a boiling stone (or a magnetic spin vane) to the dry ester. Clamping the glassware, assemble a distillation apparatus using a Hickman still and a water-cooled condenser on top of a hot plate with an aluminum heating block (Technique 14, Figure 14.5). In order to control odors, rather than to keep the reaction dry, top the apparatus with a drying tube packed loosely with a small amount of calcium chloride held in place by bits of cotton or glass wool. Begin the distillation by turning on the hot plate (about 180°C). Stir the mixture if you are using a spin vane instead of a boiling stove. Continue the distillation until only one or two drops of liquid remain in the distilling vial. If the Hickman head fills before the distillation is complete, it may be necessary to empty it using a Pasteur pipette (see Technique 14, Figure 14.6A) and transfer the distillate to a tared (preweighed) conical vial. Unless you have a sideported Hickman still, it will be necessary to remove the condenser in order to perform the transfer. When the distillation is complete, transfer the final portion of the distillate to this same vial.

Determination of Yield

Weigh the product and calculate the percentage yield of the ester. Determine its boiling point (bp 142°C) using a microscale boiling-point determination (Technique 13, Section 13.2). See the end of Experiment 14B for a spectral analysis.

14B EXPERIMENT 14B

Isopentyl Acetate (Semimicroscale Procedure)

PROCEDURE

Apparatus

Assemble a reflux apparatus on top of your hot plate using a 20- or 25-mL round-bottom flask and a water-cooled condenser (refer to Technique 7, Figure 7.6A but use a round-bottom flask instead of the conical vial). To control vapors, place a drying tube packed with calcium chloride on top of the condenser. Use a hot plate and the aluminum block with the larger set of holes for heating.

Reaction Mixture

Weigh (tare) an empty 10-mL graduated cylinder and record its weight. Place approximately 2.5 mL of isopentyl alcohol in the graduated cylinder and reweigh it to determine the weight of the alcohol. Disconnect the round-bottom flask from the reflux apparatus and transfer the alcohol into it. Do not clean or wash the graduated cylinder. Using the same graduated cylinder, measure approximately 3.5 mL of glacial acetic acid ($MW = 60$, $d = 1.06$ g/mL) and add it to the alcohol already in the flask. Using a calibrated Pasteur pipette, add 0.5 mL of concentrated sulfuric acid, mixing *immediately* (swirl), to the reaction mixture contained in the flask. Add a corundum (black) boiling stone or stirring bar and reconnect the flask. Do not use a calcium carbonate (white, marble) boiling stone, because it will dissolve in the acidic medium.

Reflux


Start water circulating in the condenser and bring the mixture to a boil. Continue heating under reflux for at least 60 minutes. Be sure to stir the mixture if you are using a stirring bar instead of a boiling stone. When the reflux period is complete, disconnect or remove the heating source and let the mixture cool to room temperature.

Extractions

Disassemble the apparatus and transfer the reaction mixture to a 15-mL capped centrifuge tube. Avoid transferring the boiling stone or stirring bar. Add 5 mL of water, cap the centrifuge tube, and mix the phases by careful shaking and venting. Allow the phases to separate and then open the cap and remove the lower aqueous layer (see a similar procedure for a conical vial in Workup). Next, extract the organic layer with 2.5 mL of aqueous sodium bicarbonate, just as you did previously with water. Extract the organic layer once again, this time with 2.5 mL of saturated aqueous sodium chloride.

Drying

Transfer the crude ester to a clean, dry, 25-mL Erlenmeyer flask and add approximately 0.5 g of anhydrous sodium sulfate. Cork the mixture and let it stand for about 10 minutes while you prepare the apparatus for distillation. If the mixture does not appear dry (the drying agent clumps and does not “flow,” the solution is cloudy, or drops of water are obvious), transfer the ester to a new, clean, dry, 25-mL Erlenmeyer flask and add a new 0.25-g portion of anhydrous sodium sulfate to complete the drying.

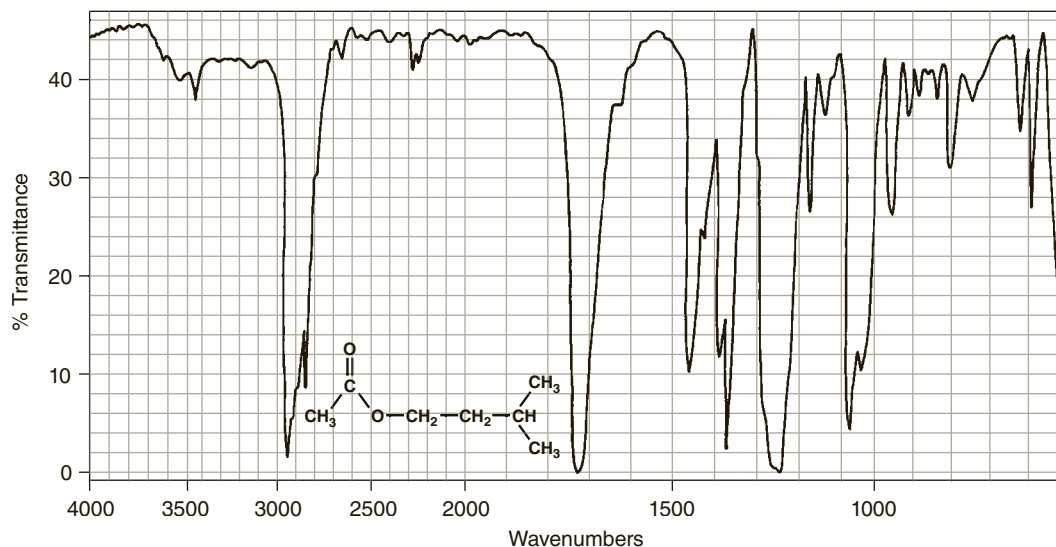

Infrared spectrum of isopentyl acetate (neat).

Distillation

Assemble a distillation apparatus using your smallest round-bottom flask to distill from (Technique 14, Figure 14.10, but insert a water condenser as shown in Experiment 7A). Use a hot plate with an aluminum block to heat. Preweigh (tare) and use a 5-mL conical vial to collect the product. Immerse the collection vial in a beaker of ice to ensure condensation and to reduce odors. Distill your ester and record its boiling-point *range* in your notebook.

Yield Determination

Weigh the product and calculate the percentage yield of the ester. At the option of your instructor, determine the boiling point using one of the methods described in Technique 13, Section 13.2. See below for a spectral analysis.

Infrared Spectroscopy

At your instructor's option, obtain an infrared spectrum using salt plates (Technique 25, Section 25.2). Compare the spectrum with the one reproduced in this experiment and include it with your report to the instructor. If any of your sample remains after performing the determination of the infrared spectrum, submit it in a properly labeled vial along with your report.

QUESTIONS

1. One method for favoring the formation of an ester is to add excess acetic acid. Suggest another method, involving the right-hand side of the equation, that will favor the formation of the ester.
2. Why is it easier to remove excess acetic acid from the products than excess isopentyl alcohol?
3. Why is the reaction mixture extracted with sodium bicarbonate? Give an equation and explain its relevance.
4. Which starting material is the limiting reagent in this procedure? Which reagent is used in excess? How great is the molar excess (how many times greater)?

5. How many grams are there in 1.00 mL of isopentyl acetate? You will need to look up the density of isopentyl acetate in a handbook.
6. How many moles of isopentyl acetate are there in 1.00 g of isopentyl acetate? You will need to calculate the molecular weight of isopentyl acetate.
7. Suppose that 1.00 mL of isopentyl alcohol was reacted with excess acetic acid and that 1.00 g of isopentyl acetate was obtained as product. Calculate the percentage yield.
8. Outline a separation scheme for isolating pure isopentyl acetate from the reaction mixture.
9. Interpret the principal absorption bands in the infrared spectrum of isopentyl acetate. (Technique 25 may be of some help in answering this question.)
10. Write a mechanism for the acid-catalyzed esterification of acetic acid with isopentyl alcohol. You may need to consult the chapter on carboxylic acids in your lecture textbook.

Terpenes and Phenylpropanoids

Anyone who has walked through a pine or cedar forest, or anyone who loves flowers and spices, knows that many plants and trees have distinctively pleasant odors. The essences or aromas of plants are due to volatile or **essential oils**, many of which have been valued since antiquity for their characteristic odors (frankincense and myrrh, for example). A list of the commercially important essential oils would run to more than 200 entries. Allspice, almond, anise, basil, bayberry, caraway, cinnamon, clove, cumin, dill, eucalyptus, garlic, jasmine, juniper, orange, peppermint, rose, sandalwood, saffron, spearmint, thyme, violet, and wintergreen are some of the familiar examples of such valuable essential oils. Essential oils are used for their pleasant odors in perfumes and incense. They are also used for their taste appeal as spices and flavoring agents in foods. A few are valued for antibacterial and antifungal action. Some are used medicinally (camphor and eucalyptus) and others as insect repellents (citronella). Chaulmoogra oil is one of the few known curative agents for leprosy. Turpentine is used as a solvent for many paint products.

Essential oil components are often found in the glands or intercellular spaces in plant tissue. They may exist in all parts of the plant but are often concentrated in the seeds or flowers. Many components of essential oils are steam-volatile and can be isolated by steam distillation. Other methods of isolating essential oils include solvent extraction and pressing (expression) methods. Esters (see the essay "Esters-Flavors and Fragrances") are frequently responsible for the characteristic odors and flavors of fruits and flowers, but other types of substances may also be important components of odor or flavor principles. Besides the esters, the ingredients of essential oils may be complex mixtures of hydrocarbons, alcohols, and carbonyl compounds. These other components usually belong to one of the two groups of natural products called **terpenes** or **phenylpropanoids**.


Terpenes

Chemical investigations of essential oils in the 19th century found that many of the compounds responsible for the pleasant odors contained exactly 10 carbon atoms. These 10-carbon compounds came to be known as terpenes if they were hydrocarbons and as **terpenoids** if they contained oxygen and were alcohols, ketones, or aldehydes.


Eventually, it was found that there are also minor and less volatile plant constituents containing 15, 20, 30, and 40 carbon atoms. Because compounds of 10 carbons were originally called terpenes, they came to be called **monoterpenes**. The other terpenes were classified in the following way.

Class	No. of Carbons	Class	No. of Carbons
Hemiterpenes	5	Diterpenes	20
Monoterpenes	10	Triterpenes	30
Sesquiterpenes	15	Tetraterpenes	40


Further chemical investigations of the terpenes, all of which contain multiples of five carbons, showed them to have a repeating structural unit based on a five-carbon pattern. This structural pattern corresponds to the arrangement of atoms in the simple five-carbon compound isoprene. Isoprene was first obtained by the thermal “cracking” of natural rubber.

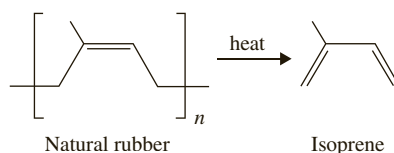

As a result of this structural similarity, a diagnostic rule for terpenes, called the **isoprene rule**, was formulated. This rule states that a terpene should be divisible, at least formally, into **isoprene units**. The structures of a number of terpenes, along with a diagrammatic division of their structures into isoprene units, is shown in the figure on the next page that accompanies this essay. Many of these compounds represent odors or flavors that should be familiar to you.

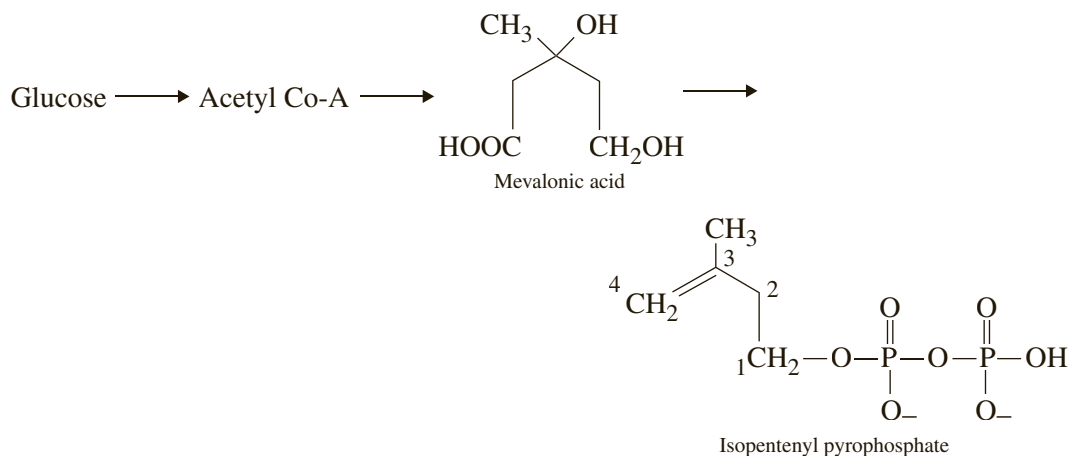


Modern research has shown that terpenes do not arise from isoprene; it has never been detected as a natural product. Instead, the terpenes arise from an important biochemical precursor compound called **mevalonic acid** (see above). This compound arises from acetyl coenzyme A, a product of the biological degradation of glucose (glycolysis), and is converted to a compound called isopentenyl pyrophosphate. Isopentenyl pyrophosphate and its isomer 3,3-dimethylallyl pyrophosphate (double bond moved to the second position) are the five-carbon building blocks used by nature to construct all the terpene compounds.

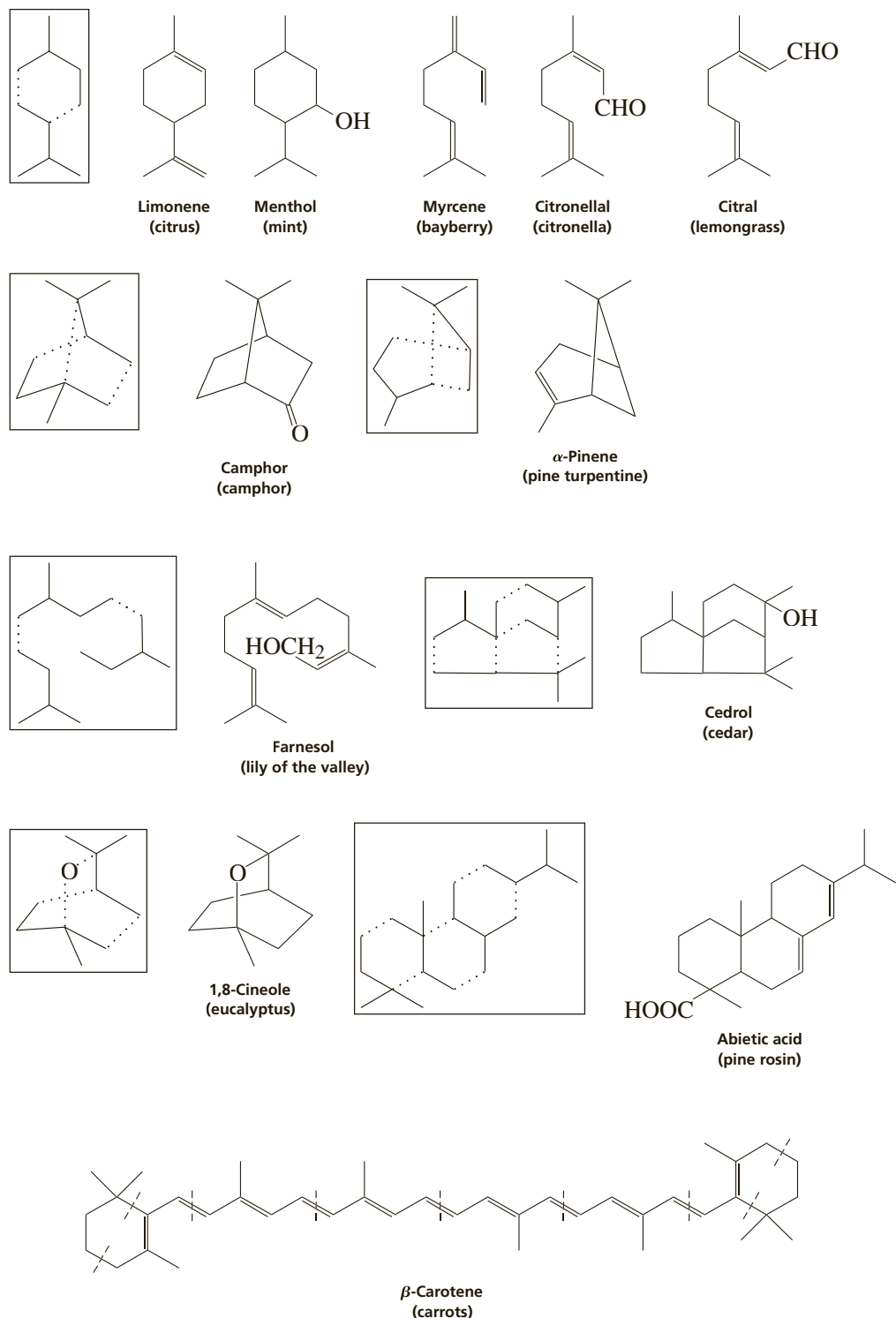


Selected terpenes.

Phenylpropanoids

Aromatic compounds, those containing a benzene ring, are also a major type of compound found in essential oils. Some of these compounds, like *p*-cymene, are actually cyclic terpenes that have been aromatized (had their ring converted to a benzene ring), but most are of a different origin.

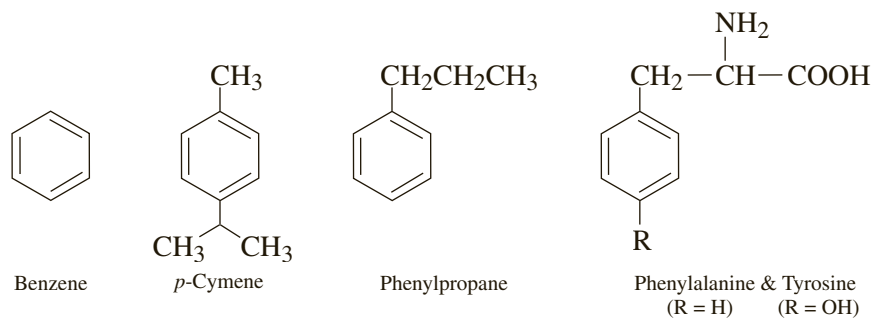

Many of these aromatic compounds are **phenylpropanoids**, compounds based on a phenylpropane skeleton. Phenylpropanoids are related in structure to the common amino acids phenylalanine and tyrosine, and many are derived from a biochemical pathway called the **shikimic acid pathway**.

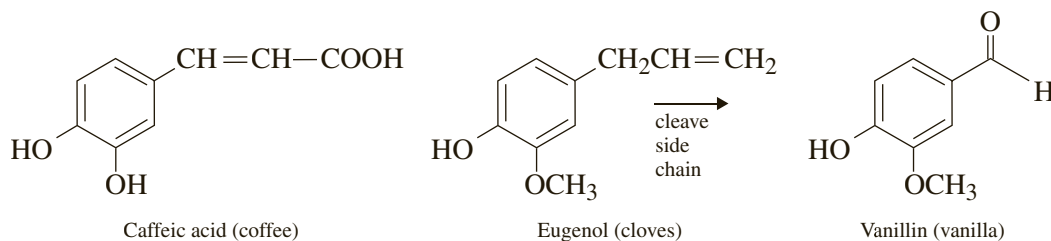

It is also common to find compounds of phenylpropanoid origin that have had the three-carbon side chain cleaved. As a result, phenylmethane derivatives, such as vanillin, are also common in plants.

REFERENCES

- Cane, D. E., Ed. Isoprenoids Including Carotenoids and Steroids. In *Comprehensive Natural Products Chemistry*, Vol. 2; Barton, D., Nakanishi, K., and Meth-Cohn, O., Eds.; Elsevier: New York, 1999.
- Cornforth, J. W. Terpene Biosynthesis. *Chem. Br.* **1968**, 4, 102.
- Geissman, T. A.; and Crout, D. H. G. *Organic Chemistry of Secondary Plant Metabolism*. Freeman, Cooper and Co.: San Francisco, 1969.
- Hendrickson, J. B. *The Molecules of Nature*. W. A. Benjamin: New York, 1965.
- Leeper, F. J.; and Vederas, J. C., Eds. *Biosynthesis: Aromatic Polyketides, Isoprenoids, Alkaloids*. Springer: New York, 2000.
- Newman, A. A. *Chemistry of Terpenes and Terpenoids: A Survey for Advanced Students and Research Workers*. Academic Press: New York, 1972.
- Pinder, A. R. *The Terpenes*. John Wiley & Sons: New York, 1960.
- Ruzicka, L. History of the Isoprene Rule. *Proc. Chem. Soc. (London)* **1959**, 341.
- Sterret, F. S. The Nature of Essential Oils, Part I. Production. *J. Chem. Edu.* **1962**, 39, 203.
- Sterret, F. S. The Nature of Essential Oils, Part II. Chemical Constituents. Analysis. *J. Chem. Edu.* **1962**, 39, 246.

Essential Oils: Extraction of Oil of Cloves by Steam Distillation


Isolation of a natural product

Steam distillation


Green chemistry

Essential oils are volatile compounds responsible for the aromas commonly associated with many plants (see essay "Terpenes and Phenylpropanoids"). The chief constituent of the essential oil from cloves is aromatic and volatile with steam. In this experiment, you will isolate the main component derived from this spice by steam distillation. Steam distillation provides a means of isolating natural products, such as essential oils, without the risk of decomposing them thermally. Identification and characterization of this essential oil will be accomplished by infrared spectroscopy.

Oil of cloves (from *Eugenia caryophyllata*) is rich in **eugenol** (4-allyl-2-methoxyphenol). Caryophyllene is present in small amounts, along with other terpenes. Eugenol (bp 250°C) is a phenol, or an aromatic hydroxy compound.

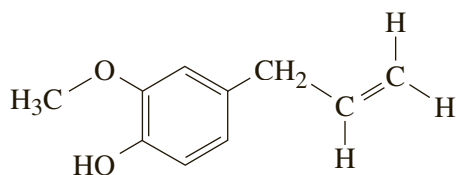

Eugenol

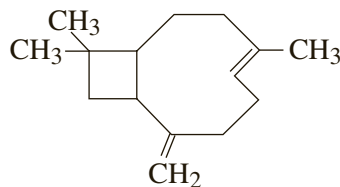

Caryophyllene

REQUIRED READING

Review: Techniques 5 and 6

Technique 7 Reaction Methods, Section 7.10

Technique 12 Extractions, Separations, and Drying Agents, Sections 12.4 and 12.9

Technique 25 Infrared Spectroscopy

New: Technique 18 Steam Distillation

Essay Terpenes and Phenylpropanoids

If performing the optional proton NMR analysis, also read:

Technique 26 Nuclear Magnetic Resonance Spectroscopy

SPECIAL INSTRUCTIONS

Be careful when handling methylene chloride. It is a toxic solvent, and you should not breathe it excessively or spill it on yourself. To complete the distillation in a reasonable time, boil the mixture as rapidly as possible without allowing the boiling mixture to rise above the neck of the Hickman head. This requires that you work with careful attention during the distillation procedure. The distillation requires 1–2 hours.

SUGGESTED WASTE DISPOSAL

You must dispose of methylene chloride in a waste container marked for the disposal of halogenated organic waste. Any residue from the ground cloves can be disposed of in an ordinary trash can. Any aqueous solutions should be placed in the container specially designated for aqueous wastes.

15A

 EXPERIMENT 15A

Oil of Cloves (Microscale Procedure)

PROCEDURE

Assemble a steam distillation apparatus, as shown in Technique 18, Figure 18.3. Be sure to include the water condenser, as shown in the illustration. Use a 20- or 25-mL round-bottom flask as a distillation flask and either an aluminum block or a sand bath to heat the distillation flask. If you use a sand bath, you may need to cover the sand bath and distillation flask with aluminum foil.

Weigh approximately 1.0 g of ground cloves or clove buds onto a weighing paper and record the exact weight. If your spice is already ground, you may proceed without grinding it; if you use clove buds, cut the buds into small pieces. Mix the spice with 12–15 mL of water in the round-bottom distillation flask, add a magnetic stirring bar, and attach the flask to the distillation apparatus. Allow the spice to soak in the water for about 15 minutes before beginning the heating. Be sure that all the spice is thoroughly wetted. Swirl the flask gently, if needed.

Steam Distillation

Turn on the cooling water in the condenser, begin stirring the mixture in the distillation flask, and begin heating the mixture to provide a steady rate of distillation. The temperature for the heating device should be about 130°C. If you approach the boiling point too quickly, you may have difficulty with frothing or bump-over. You need to find the amount of heating that provides a steady rate of distillation but avoids frothing or bumping. A good rate of distillation would be to have one drop of liquid collected every 2–5 seconds.

As the distillation proceeds, use a Pasteur pipette (5¾-inch) to transfer the distillate from the reservoir of the Hickman head to a 15-mL screw-cap centrifuge tube. If you are using a Hickman head with a side port, you can easily remove the distillate by opening the side port and withdrawing the liquid. If your Hickman head

does not have a side port, you need to remove the condenser from the top of the distillation apparatus to remove the distillate. In this case, the transferring operation is best accomplished if the Pasteur pipette is bent slightly at the end. Continue distillation until 5–8 mL of distillate has been collected.

Normally in a steam distillation, the distillate is somewhat cloudy owing to separation of the essential oil as the vapors cool. You may not notice this, but you will still obtain satisfactory results.

Extraction of the Essential Oil

Collect all the distillate in a 15-mL screw-cap centrifuge tube. Using a calibrated Pasteur pipette (see Experiment 1), add 2.0 mL of methylene chloride (dichloromethane) to extract the distillate. Cap the tube securely and shake it vigorously with frequent venting. Allow the layers to separate. Using a Pasteur pipette, transfer the lower methylene chloride layer to a clean, dry, 5-mL conical vial. Repeat this extraction procedure two more times with fresh 1.0-mL portions of methylene chloride and combine all the methylene chloride extracts in the same 5-mL conical vial that you used for the first extraction. If there are drops of water in the vial, it will be necessary to transfer the methylene chloride solution with a dry Pasteur pipette to another dry conical vial.

Drying

Dry the methylene chloride solution by adding granular anhydrous sodium sulfate to the conical vial (see Technique 12, Section 12.9). Let the solution stand for 10–15 minutes with occasional stirring.

Evaporation

While the organic solution is being dried, clean and dry a 5-mL conical vial and weigh (tare) it accurately. With a clean, dry Pasteur pipette, transfer the dried organic layer to the tared vial, leaving the drying agent behind. Use small amounts of methylene chloride to rinse the solution completely into the tared vial. Be careful to keep any of the sodium sulfate from being transferred. Working in a hood, evaporate the methylene chloride from the solution by using a gentle stream of dry air or nitrogen while heating the vial in a warm water bath (temperature about 40°C). (See Technique 7, Section 7.10.) It is important that the stream of air or nitrogen be gentle or you will force your solution out of the conical vial. In addition, be careful not to overheat the sample. Be careful not to continue the evaporation beyond the point where all the methylene chloride has evaporated. Your product is a volatile oil (i.e., a liquid), and if you continue to heat and evaporate the liquid beyond the point where the solvent has been removed, you will likely lose your sample.

Yield Determination

When the solvent has been removed, weigh the conical vial. Calculate the weight percentage recovery (see Technique 2, Section 2.2C) of the oil from the original amount of spice used. See Section 15B for a spectral analysis.

15B EXPERIMENT 15B*Oil of Cloves (Semimicroscale Procedure)***PROCEDURE***Apparatus*

Assemble a semimicroscale distillation apparatus, as shown in Technique 14, Figure 14.10. Use a 20- or 25-mL round-bottom flask as the distillation flask and either an aluminum block or a sand bath to heat the distillation flask. If you use a sand bath, you may need to cover the sand bath and distillation flask with aluminum foil.

Preparation


Use the amounts of cloves and water described in Experiment 15A.

Steam Distillation

Proceed with the distillation as described in Experiment 15A. Note, however, that you will not have to remove distillate during the course of the distillation. Continue with the extraction, drying, evaporation, and yield determination, as described in Experiment 15A.

**Spectroscopy
(Experiment 14A and
14B)***Infrared Spectrum*

Obtain the infrared spectrum of the oil as a pure liquid sample (Technique 25, Section 25.2). Small amounts of water will damage the salt plates that are used as cells in infrared spectroscopy.

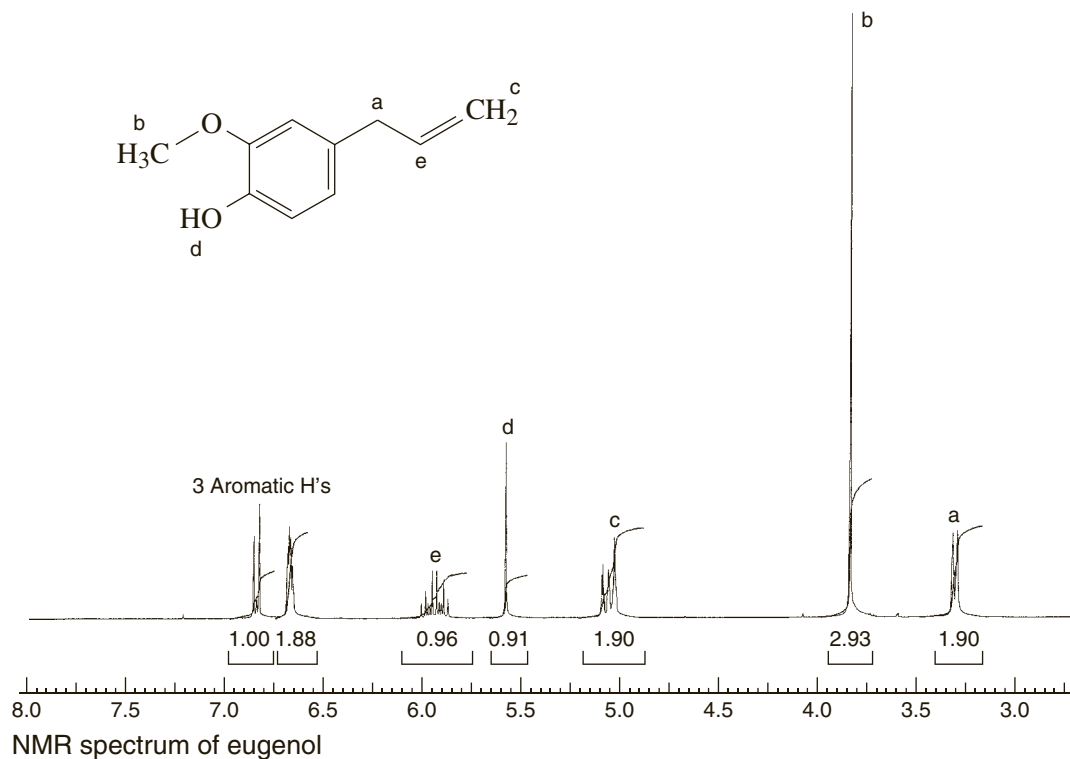

NOTE: Before proceeding with infrared spectroscopy, check with your instructor to make sure your sample is acceptable.

It may be necessary to use a capillary tube to transfer a sufficient amount of liquid to the salt plates. If the amount of liquid is too small to transfer, add one or two drops of methylene chloride to aid in the transfer. Gently blow on the plate to evaporate the solvent. Include the infrared spectrum in your laboratory report, along with an interpretation of the principal absorption peaks.

NMR Spectrum

At the instructor's option, determine the nuclear magnetic resonance spectrum of the oil (Technique 26, Section 26.1).

QUESTIONS

1. Why is eugenol steam-distilled rather than purified by simple distillation?
2. A natural product ($MW = 150$) distills with steam at a boiling temperature of 99°C at atmospheric pressure. The vapor pressure of water at 99°C is 733 mm Hg.
 - a. Calculate the weight of the natural product that codistills with each gram of water at 99°C .
 - b. How much water must be removed by steam distillation to recover this natural product from 0.5 g of a spice that contains 10% of the desired substance?
3. In a steam distillation, the amount of water actually distilled is usually greater than the amount calculated, assuming that both water and organic substance exert the same vapor pressure when they are mixed that they exert when each is pure. Why does one recover more water in the steam distillation than was calculated? (*Hint:* Are the organic compound and water truly immiscible?)
4. Explain how caryophyllene fits the isoprene rule (see essay, "Terpenes and Phenylpropanoids").

Stereochemical Theory of Odor

The human nose has an almost unbelievable ability to distinguish odors. Just consider for a few moments the different substances you can recognize by odor alone. Your list should be long. A person with a trained nose, a perfumer, for instance, can often recognize even individual components in a mixture. Who has not met at least one cook who could sniff almost any culinary dish and identify the seasonings and spices that were used? The olfactory centers in the nose can identify odorous substances even in small amounts. Studies have shown that with some substances, as little as one 10 millionth of a gram (10^{-7} g) can be perceived. Many animals, for example, dogs and insects, have an even lower threshold of smell than humans do (see essay on pheromones that precedes Experiment 50).

There have been many theories of odor, but few have persisted. Strangely enough, one of the oldest theories, although in modern dress, is still the most current theory. Lucretius, one of the early Greek atomists, suggested that substances having odor gave off a vapor of tiny "atoms," all of the same shape and size, and that they gave rise to the perception of odor when they entered pores in the nose. The pores would have to be of various shapes, and the odor perceived would depend on which pores the atoms were able to enter. We now have many similar theories about the action of drugs (receptor-site theory) and the interaction of enzymes with their substrates (the lock-and-key hypothesis).

A substance must have certain physical characteristics to have the property of odor. First, it must be volatile enough to give off a vapor that can reach the nostrils. Second, once it reaches the nostrils, it must be somewhat water soluble, even if only to a small degree, so that it can pass through the layer of moisture (mucus) that covers the nerve endings in the olfactory area. Third, it must have lipid solubility to allow it to penetrate the lipid (fat) layers that form the surface membranes of the nerve cell endings.

Once we pass these criteria, we come to the heart of the question. Why do substances have different odors? In 1949, R. W. Moncrieff, a Scot, resurrected Lucretius' hypothesis. He proposed that in the olfactory area of the nose is a system of receptor cells of several types and shapes. He further suggested that each receptor site corresponded to a different type of primary odor. Molecules that would fit these receptor sites would display the characteristics of that primary odor. It would not be necessary for the entire molecule to fit into the receptor, so for larger molecules, any portion might fit into the receptor and activate it. Molecules having complex odors would presumably be able to activate several types of receptors.

Moncrieff's hypothesis was strengthened substantially by the work of J. E. Amoore, who began studying the subject as an undergraduate at Oxford in 1952. After an extensive search of the chemical literature, Amoore concluded that there were only seven basic primary odors. By sorting molecules with similar odor types, he even formulated possible shapes for the seven necessary receptors. For instance, from the literature he culled more than 100 compounds that were described as having

a “camphoraceous” odor. Comparing the sizes and shapes of all these molecules, he postulated a three-dimensional shape for a camphoraceous receptor site. Similarly, he derived shapes for the other six receptor sites. The seven primary receptor sites he formulated are shown in Figure 1, along with a typical prototype molecule of the appropriate shape to fit the receptor. The shapes of the sites are shown in perspective. Pungent and putrid odors were not thought to require a particular shape in the odorless molecules but rather to need a particular type of charge distribution.

You can verify quickly that compounds with molecules of roughly similar shape have similar odors if you compare nitrobenzene and acetophenone with benzaldehyde or *d*-camphor and hexachloroethane with cyclooctane. Each group of substances has the same basic odor *type* (primary), but the individual molecules differ in the *quality* of the odor. Some of the odors are sharp, some pungent, others sweet, and so on. The second group of substances all have a camphoraceous odor, and the molecules of these substances all have approximately the same shape.

An interesting corollary to the Amoore theory is the postulate that if the receptor sites are chiral, then optical isomers (enantiomers) of a given substance might have *different* odors. This circumstance proves true in several cases. It is true for (+)- and (–)-carvone; we investigate the idea in Experiment 16 in this textbook.

The theory changed dramatically in 1991 because of the biochemical research of Richard Axel and Linda Buck, who was a postdoctoral student in Axel’s research group. Subsequently, Buck founded her own group that also continued research on the nature of the sense of smell. In 2004, Axel and Buck won the Nobel Prize in Physiology or Medicine for their combined work during the previous decade.

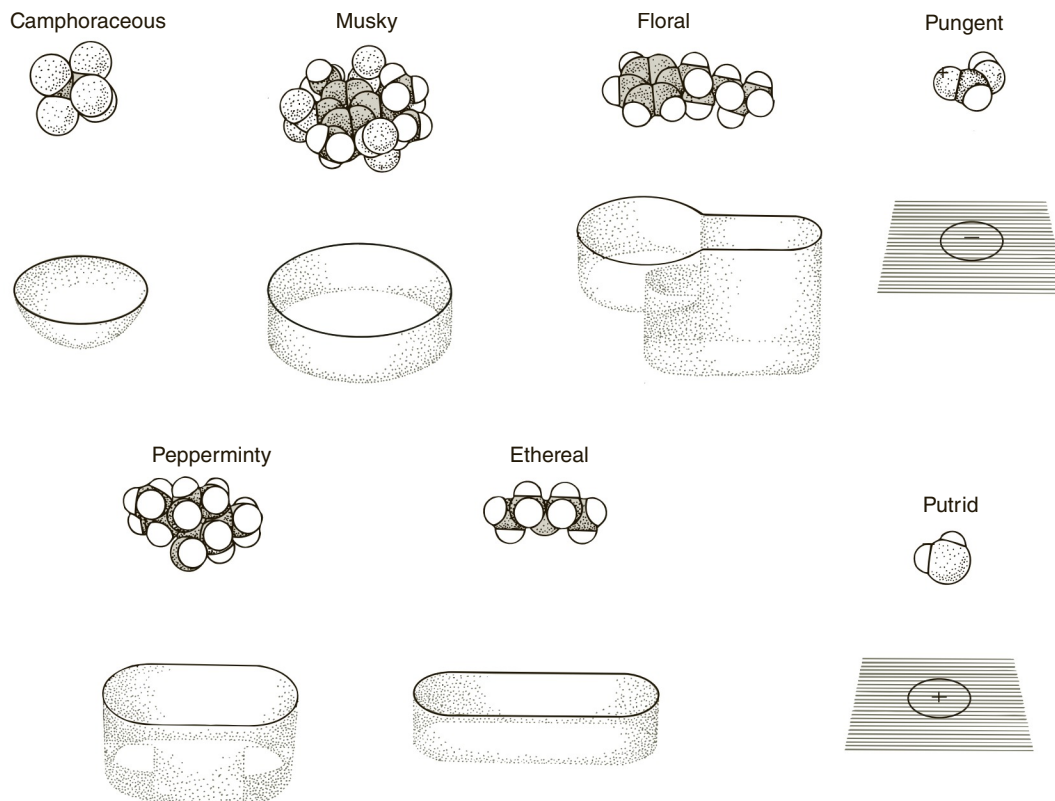


Figure 1
Seven Primary Odor Receptor Sites.

From “The Stereochemical Theory of Odor,” by J. E. Amoore, J. W. Johnston, Jr., and M. Rubin. *Scientific American*, 210:42–49 Copyright © 1964 by Scientific American, Inc. All rights reserved. Reprinted by permission.

The 1991 paper, working with mice, described a family of membrane-spanning receptor proteins found in a small area of the upper nose called the olfactory epithelium. Mice have genes that can encode as many as 1000 types of receptor proteins. Subsequent work has estimated that humans, who have a lesser-developed sense of smell than mice, encode only about 350 of these receptor proteins. Each of these protein receptors is located on the surface of the olfactory epithelium and is connected to a single nerve cell (neuron) located in the epithelium. The neuron “fires” or sends a signal when an odorant molecule binds to the active site of the protein. The signal is carried across the bones of the skull and into a node in an area of the brain called the olfactory bulb. The signals from all receptors are processed in the olfactory bulb and sent to the memory area of the brain where recognition of the odor takes place. Figure 2 shows a schematic of the olfactory region.

The signals from all of the types of protein receptors are collected, or integrated, in the olfactory bulb. The node (a postulated feature) is a common connection where the signals from each type of cell are collected and sent to memory, each with an intensity proportional to the numbers of cells that were stimulated by the odorant molecules. Because a given odorant molecule should be capable of binding to more than one type of receptor and because many odors are composed of more than one type of molecule, the signal sent to memory should be a complex combinatorial pattern consisting of contributions from several nodes, each with a different intensity value. This system should allow a human to recognize as many as 10,000 odors and for mice to recognize many more. The memory region in the brain can also make associations based on a given pattern. For instance, cinnamaldehyde can be recognized as the odor of the spice cinnamon, but it can also be associated with other items such as apple pie, cinnamon rolls, apple strudel, spiced cider, and, of course, pleasure. A figure showing these associations, but with a limitation of only a few receptors represented, is shown in Figure 3.

Although our modern understanding of the detection of odor has evolved to become a more highly detailed theory than the one proposed by Lucretius, it would appear that his fundamental hypothesis was correct and has even withstood the scrutiny of modern science.

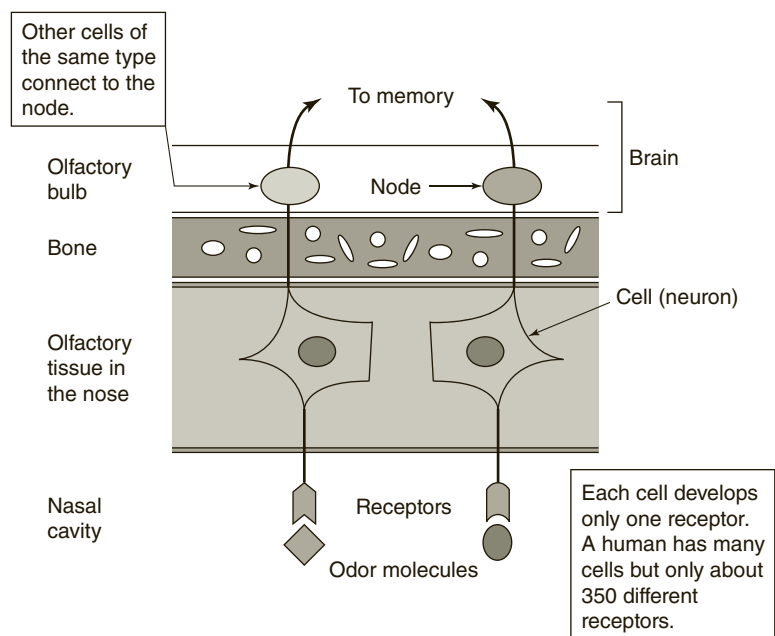


Figure 2
Odor receptors in the nose.

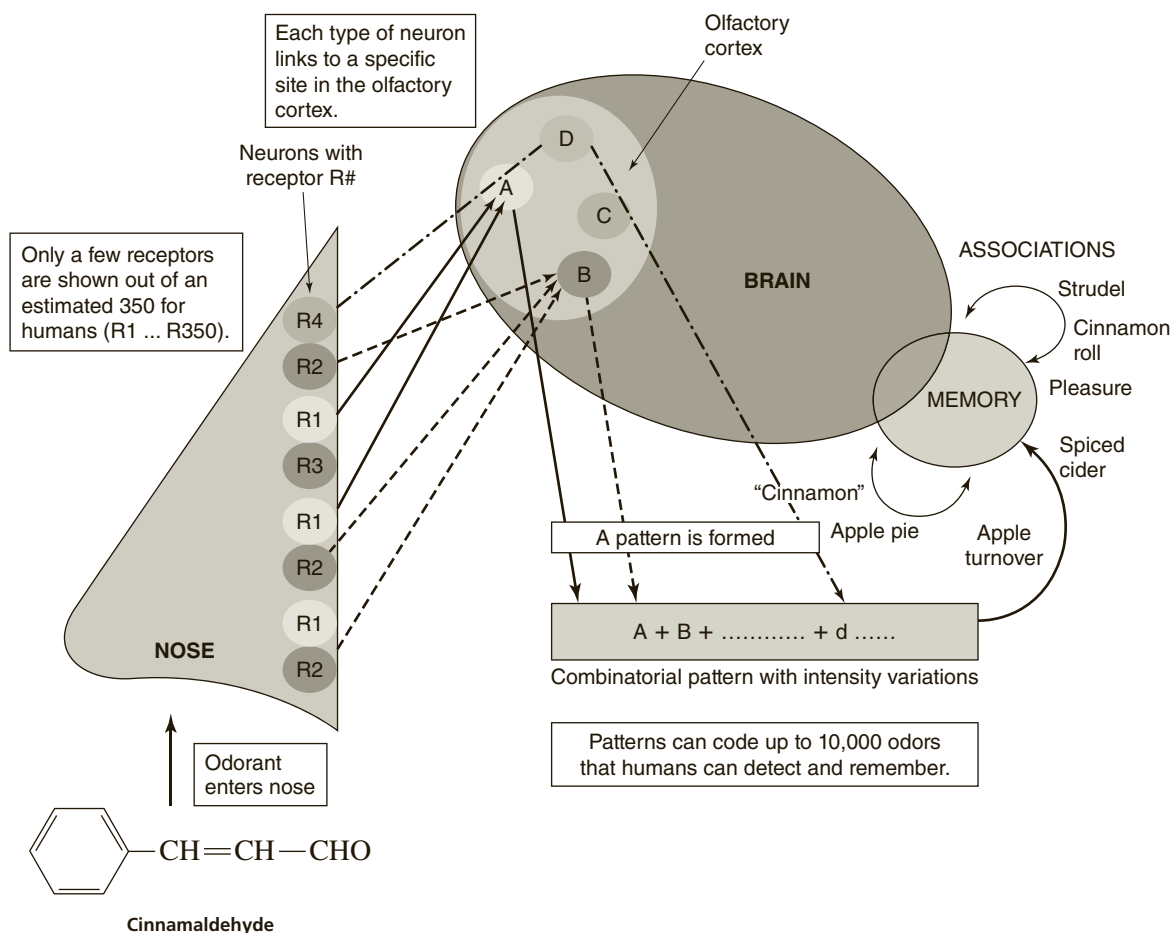


Figure 3
Nobel prize theory of the detection of odors.

REFERENCES

- Amoore, J. E.; Johnson, J. W., Jr.; and Rubin, M. The Stereochemical Theory of Odor. *Sci. Am.* **1964**, 210 (Feb), 1.
- Amoore, J. E.; Johnson, J. W., Jr.; and Rubin, M. The Stereochemical Theory of Olfaction. *Proc. Sci. Sec. TGA* (Special Supplement to No. 37) **1962** (Oct), 1–47.
- Buck, L. The Molecular Architecture of Odor and Pheromone Sensing in Mammals. *Cell* **2000**, 100(6) (Mar), 611–618.
- Buck, L.; and Axel, R. A Novel Multigene Family May Encode Odorant Receptors: A Molecular Basis for Odor Recognition. *Cell* **1991**, 65(1) (Apr): 175–187.
- Lipkowitz, K. B. Molecular Modeling in Organic Chemistry: Correlating Odors with Molecular Structure. *J. Chem. Edu.* **1989**, 66 (Apr), 275.
- Malnic, B.; Hirono, J.; Sato, T.; and Buck, L. Combinatorial Receptor Codes for Odors. *Cell* **1999**, 96(5) (Mar), 713–723.
- Moncrieff, R. W. *The Chemical Senses*. Routledge & Kegan Paul: London, 1976.
- Roderick, W. R. Current Ideas on the Chemical Basis of Olfaction. *J. Chem. Edu.* **1966**, 43 (Oct), 510–519.
- Zou, Z.; Horowitz, L.; Montmayeur, J.; Snapper, S.; and Buck, L. Genetic Tracing Reveals a Stereotyped Sensory Map in the Olfactory Cortex. *Nature* **2001**, 414(6843) (Nov), 173–179.

Spearmint and Caraway Oil: (+)- and (-)-Carvones


Stereochemistry

Gas chromatography

Polarimetry

Spectroscopy


Refractometry

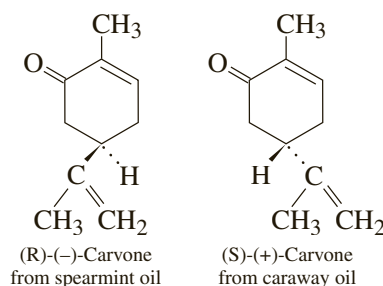


In this experiment, you will compare (+)-carvone from caraway oil to (-)-carvone from spearmint oil, using gas chromatography. If you have the proper preparative-scale gas-chromatographic equipment, it should be possible to prepare pure samples of each of the carvones from their respective oils. If this equipment is not available, the instructor will provide pure samples of the two carvones obtained from a commercial source, and any gas-chromatographic work will be strictly analytical.

The odors of the two enantiomeric carvones are distinctly different from each other. The presence of one or the other isomer is responsible for the characteristic odors of each oil. The difference in the odors is to be expected because the odor receptors in the nose are chiral (see essay, "Stereochemical Theory of Odor"). This phenomenon, in which a chiral receptor interacts differently with each enantiomer of a chiral compound, is called **chiral recognition**.

Although we should expect the optical rotations of the isomers (enantiomers) to be of opposite sign, the other physical properties should be identical. Thus, for both (+)- and (-)-carvone, we predict that the infrared and nuclear magnetic resonance spectra, the gas-chromatographic retention times, the refractive indices, and the boiling points will be identical. Hence, the only differences in properties you will observe for the two carvones are the odors and the signs of rotation in a polarimeter.

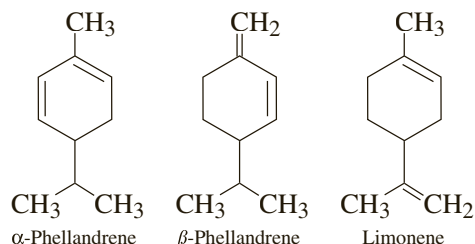

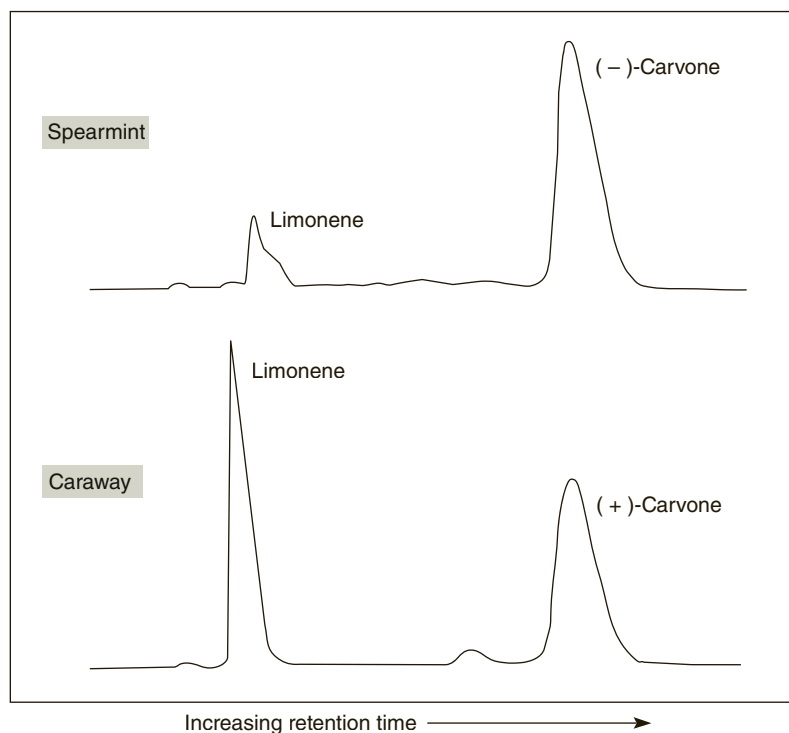

Gas chromatograms of caraway and spearmint oil.

Caraway oil contains mainly limonene and (+)-carvone. The gas chromatogram for this oil is shown in the figure. The (+)-carvone (bp 203°C) can easily be separated from the lower-boiling limonene (bp 177°C) by gas chromatography, as shown in the figure. If one has a preparative gas chromatograph, the (+)-carvone and limonene can be collected separately as they elute from the gas chromatography column. **Spearmint oil** contains mainly (–)-carvone with a smaller amount of limonene and very small amounts of the lower-boiling terpenes, α - and β -phellandrene. The gas chromatogram for this oil is also shown in the figure. With preparative equipment, you can easily collect the (–)-carvone as it exits the column. It is more difficult, however, to collect limonene in a pure form. It is likely to be contaminated with the other terpenes because they all have similar boiling points.

REQUIRED READING

<i>Review:</i>	Experiment 1	Introduction to Microscale Laboratory
	Technique 25	Infrared Spectroscopy
<i>New:</i>	Technique 22	Gas Chromatography
	Technique 23	Polarimetry
	Essay	Stereochemical Theory of Odor

If performing any of the optional procedures, read as appropriate:

Technique 13	Physical Constants of Liquids, Boiling Points
Technique 24	Refractometry
Technique 26	Nuclear Magnetic Resonance Spectroscopy
Technique 27	Carbon-13 Nuclear Magnetic Resonance Spectroscopy

SPECIAL INSTRUCTIONS

Your instructor will either assign you spearmint or caraway oil or have you choose one. You will also be given instructions on which procedures from Part A you are to perform. You should compare your data with those of someone who has studied the other enantiomer.

NOTE: If a gas chromatograph is not available, this experiment can be performed with spearmint and caraway oils and pure commercial samples of the (+)- and (-)-carvones.

If the proper equipment is available, your instructor may require you to perform a gas-chromatographic analysis. If preparative gas chromatography is available, you will be asked to isolate the carvone from your oil (Part B). Otherwise, if you are using analytical equipment, you will be able to compare only the retention times and integrals from your oil to those of the other essential oil.

Although preparative gas chromatography will yield enough sample to do spectra, it will not yield enough material to do the polarimetry. Therefore, if you are required to determine the optical rotation of the pure samples, whether or not you perform preparative gas chromatography, your instructor will provide a pre-filled polarimeter tube for each sample.

NOTES TO THE INSTRUCTOR

This experiment may be scheduled along with another experiment. It is best if students work in pairs, each student using a different oil. An appointment schedule for using the gas chromatograph should be arranged so that students are able to make efficient use of their time. You should prepare chromatograms using both carvone isomers and limonene as reference standards. Appropriate reference standards include a mixture of (+)-carvone and limonene and a second mixture of (-)-carvone and limonene. The chromatograms should be posted with retention times, or each student should be provided with a copy of the appropriate chromatogram.

The gas chromatograph should be prepared as follows: column temperature, 200°C; injection and detector temperature, 210°C; carrier gas flow rate, 20 mL/min. The recommended column is 8 feet long, with a stationary phase such as Carbowax 20M. It is convenient to use a Gow-Mac 69-350 instrument with the preparative accessory system for this experiment.

You should fill polarimeter cells (0.5 dm) in advance with the undiluted (+)- and (-)-carvones. There should also be four bottles containing spearmint and caraway oils and (+)- and (-)-carvone. Both enantiomers of carvone are commercially available.

PROCEDURE

Part A. Analysis of the Carvones

The samples (either those obtained from gas chromatography, Part B, or commercial samples) should be analyzed by the following methods. The instructor will indicate which methods to use. Compare your results with those obtained by someone who used a different oil. In addition, measure the observed rotation of the commercial samples of (+)-carvone and (-)-carvone. The instructor will supply pre-filled polarimeter tubes.

Analyses to Be Performed on Spearmint and Caraway Oils:

Odor Carefully smell the containers of spearmint and caraway oil and of the two carvones. About 8–10% of the population cannot detect the difference in the odors of the optical isomers. Most people, however, find the difference quite obvious. Record your impressions.

Analytical Gas Chromatography If you separated your sample by preparative gas chromatography in Part B, you should already have your chromatogram. In this case, you should compare it to one done by someone using the other oil. Be sure to obtain retention times and integrals or obtain a copy of the other person's chromatogram.

If you did not perform Part B, obtain the analytical gas chromatograms of your assigned oil—spearmint or caraway—and obtain the result from the other oil from someone else. The instructor may prefer to perform the sample injections or have a laboratory assistant perform them. The sample injection procedure requires careful technique, and the special microliter syringes that are required are delicate and expensive. If you are to perform the injections yourself, your instructor will give you adequate instruction beforehand.

For both oils, determine the retention times of the components (see Technique 22, Section 22.7). Calculate the percentage composition of the two essential oils by one of the methods explained in the same section.

Analyses to Be Performed on the Purified Carvones:

Polarimetry With the help of the instructor or assistant, obtain the observed optical rotation α of the pure (+)-carvone and (–)-carvone samples. These are provided in prefilled polarimeter tubes. The specific rotation $[\alpha]_D$ is calculated from the relationship given in Technique 23, Section 23.2. The concentration c will equal the density of the substances analyzed at 20°C. The values, obtained from actual commercial samples, are 0.9608 g/mL for (+)-carvone and 0.9593 g/mL for (–)-carvone. The literature values for the specific rotations are as follows: $[\alpha]_D^{20} = +61.7^\circ$ for (+)-carvone and -62.5° for (–)-carvone. These values are not identical, because trace amounts of impurities are present.

Polarimetry does not work well on the crude spearmint and caraway oils, because large amounts of limonene and other impurities are present.


Infrared Spectroscopy Obtain the infrared spectrum of the (–)-carvone sample from spearmint or of the (+)-carvone sample from caraway (see Technique 25, Section 25.2). Compare your result with that of a person working with the other isomer. At the option of the instructor, obtain the infrared spectrum of the (+)-limonene, which is found in both oils. If possible, determine all spectra using neat samples. If you isolated the samples by preparative gas chromatography, it may be necessary to add one to two drops of carbon tetrachloride to the sample. Thoroughly mix the liquids by drawing the mixture into a Pasteur pipette and expelling several times. It may be helpful to draw the end of the pipette to a narrow tip in order to withdraw all the liquid in the conical vial. As an alternative, use a microsyringe. Obtain a spectrum on this solution, as described in Technique 25, Section 25.2.

Nuclear Magnetic Resonance Spectroscopy Using an NMR instrument, obtain a proton NMR spectrum of your carvone. Compare your spectrum with the NMR spectra for (–)-carvone and (+)-limonene shown in this experiment. Attempt to assign as many peaks as you can. If your NMR instrument is capable of obtaining a carbon-13


NMR spectrum, determine a carbon-13 spectrum. Compare your spectrum of carvone with the carbon-13 NMR spectrum shown in this experiment. Once again, attempt to assign the peaks.

Boiling Point Determine the boiling point of the carvone you were assigned. Use the microboiling-point technique (Technique 13, Section 13.2). The boiling points for both carvones are 230°C at atmospheric pressure. Compare your result to that of someone using the other carvone.


Refractive Index Use the technique for obtaining the refractive index on a small volume of liquid, as described in Technique 24, Section 24.2. Obtain the refractive index for the carvone you separated (Part B) or for the one assigned. Compare your value to that obtained by someone using the other isomer. At 20°C, the (+)- and (-)-carvones have the same refractive index, equal to 1.4989.

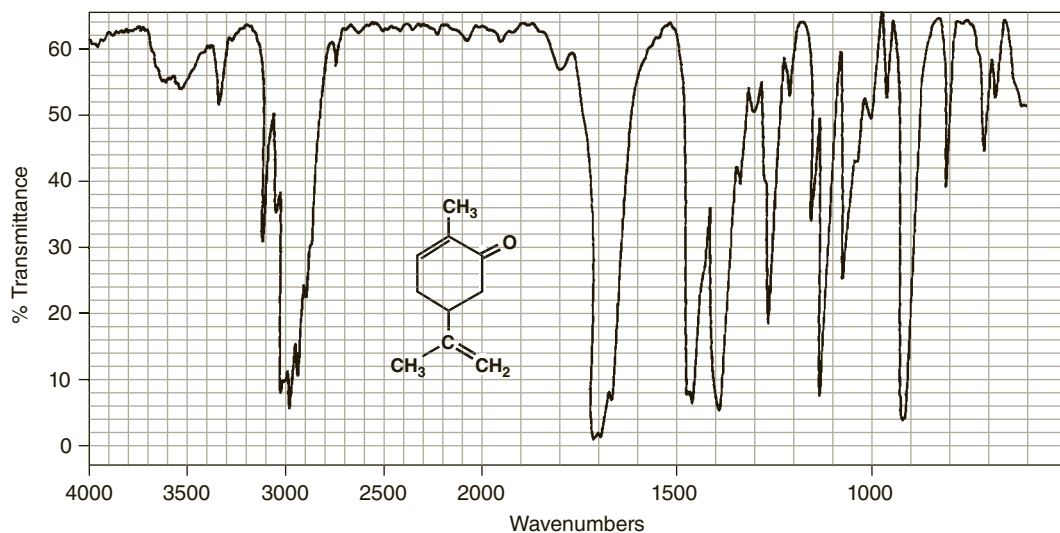


Infrared Spectrum of carvone (neat).

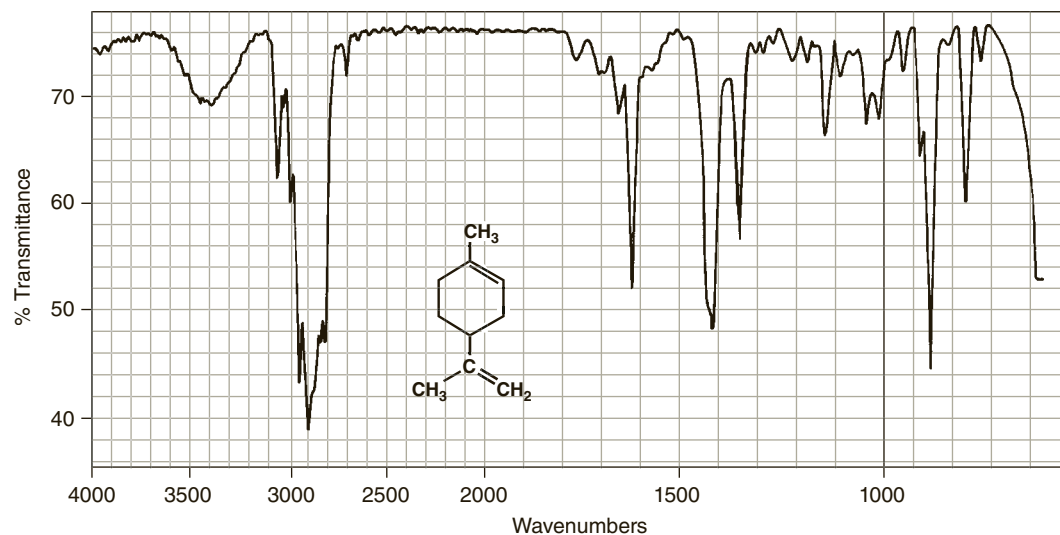


Infrared spectrum of limonene (neat).

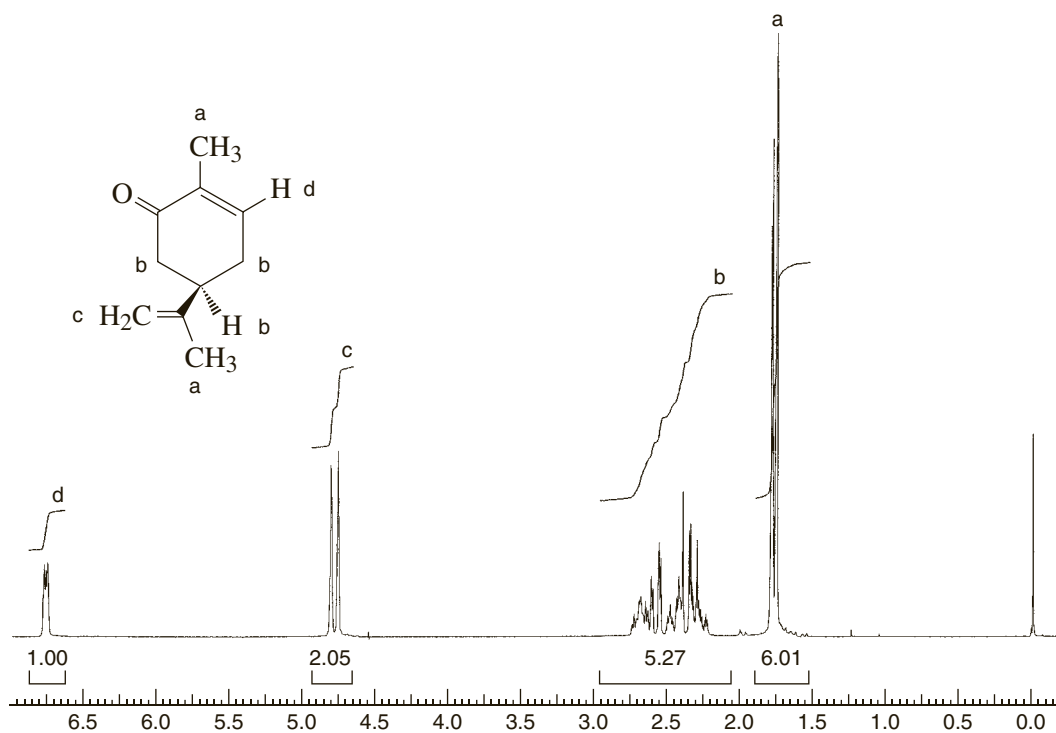

NMR spectrum of (-)-carvone from spearmint oil.

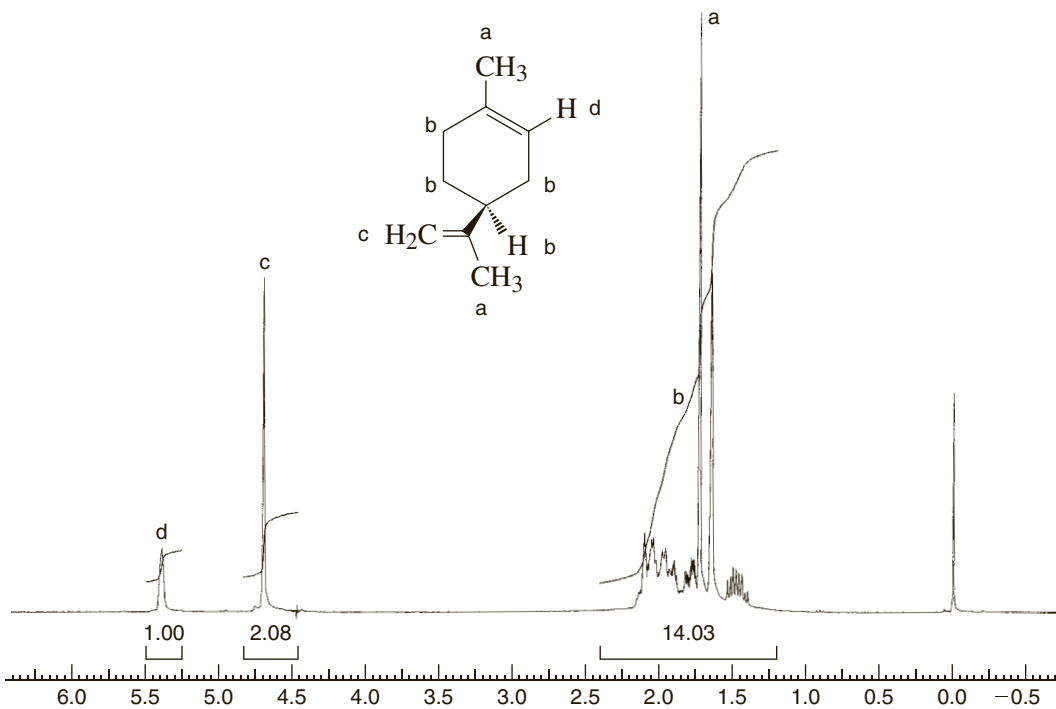

NMR spectrum of (+)-limonene.

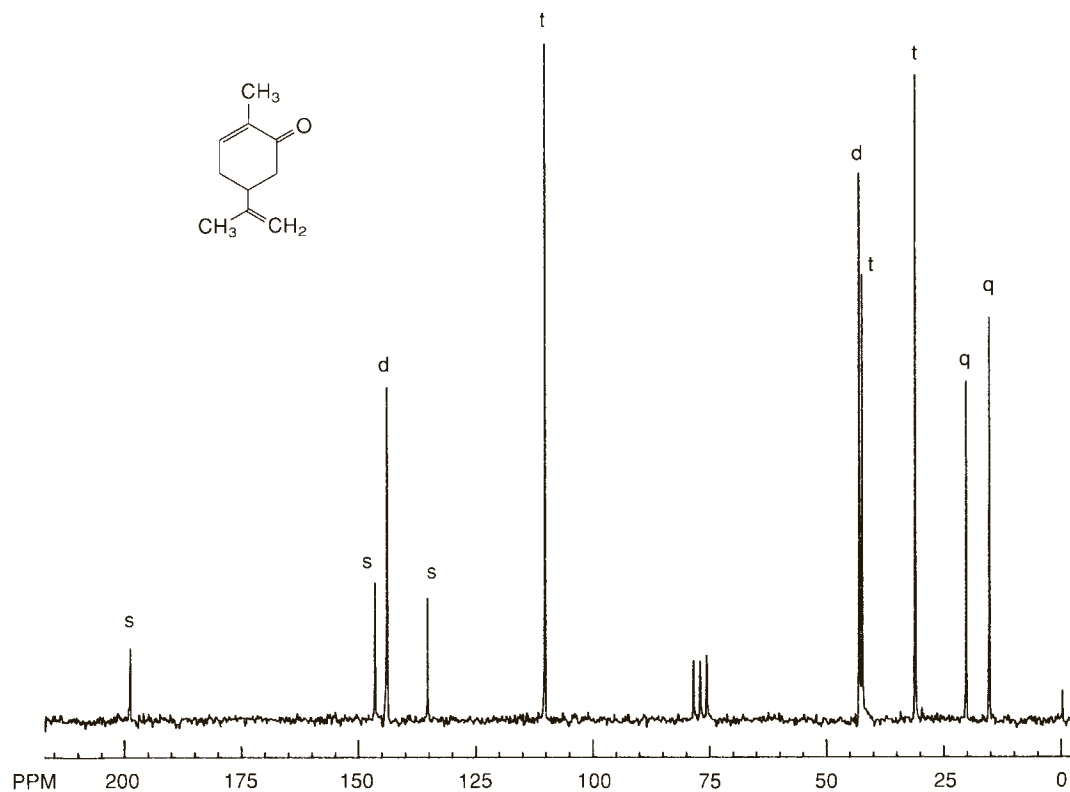

Decoupled carbon-13 spectrum of carvone, CDCl_3 . Letters indicate appearance of spectrum when carbons are coupled to protons (s = singlet, d = doublet, t = triplet, q = quartet).

Part B. Separation by Gas Chromatography (Optional)

The instructor may prefer to perform the sample injections or have a laboratory assistant perform them. The sample injection procedure requires careful technique, and the special microliter syringes that are required are delicate and expensive. If you are to perform the sample injections, your instructor will give you adequate instruction beforehand.

Inject 50 μL of caraway or spearmint oil on the gas-chromatography column. Just before a component of the oil (limonene or carvone) elutes from the column, install a gas-collection tube at the exit port, as described in Technique 22, Section 22.11. To determine when to connect the gas-collection tube, refer to the chromatograms prepared by your instructor. These chromatograms have been run on the same instrument you are using under the same conditions. Ideally, you should connect the gas-collection tube just before the limonene or carvone elutes from the column and remove the tube as soon as all the component has been collected but before any other compound begins to elute from the column. You can accomplish this most easily by watching the recorder as your sample passes through the column. The collection tube is connected (if possible) just before a peak is produced or as soon as a deflection is observed. When the pen returns to the baseline, remove the gas collection tube.

This procedure is relatively easy for collecting the carvone component of both oils and for collecting the limonene in caraway oil. Because of the presence of several terpenes in spearmint oil, it is somewhat more difficult to isolate a pure sample of limonene from spearmint oil (see the chromatogram in the introductory section of this experiment). In this case, you must try to collect only the limonene component and not any other compounds, such as the terpene, which produces a shoulder on the limonene peak in the chromatogram for spearmint oil.

After collecting the samples, insert the ground joint of the collection tube into a 0.1-mL conical vial, using an O-ring and screw cap to fasten the two pieces together securely. Place this assembly into a test tube, as shown in Technique 22, Figure 22.11. Put cotton on the bottom of the tube and use a rubber septum cap to hold the assembly in place and to prevent breakage. Balance the centrifuge by placing a tube of equal weight on the opposite side (this could be your other sample or someone else's sample). During centrifugation, the sample is forced into the bottom of the conical vial. Disassemble the apparatus, cap the vial, and perform the analyses described in Part A. You should have enough sample to perform the infrared and NMR spectroscopy, but your instructor may need to provide additional sample to perform the other procedures.

REFERENCES

- Friedman, L., and Miller, J. G. Odor, Incongruity, and Chirality. *Science*, 172 (1971): 1044.
- Murov, S. L., and Pickering, M. The Odor of Optical Isomers. *Journal of Chemical Education*, 50 (1973): 74.
- Russell, G. F., and Hills, J. I. Odor Differences Between Enantiomeric Isomers. *Science*, 172 (1971): 1043.

QUESTIONS

1. Interpret the infrared spectra for carvone and limonene and the proton and carbon-13 NMR spectra of carvone.
2. Identify the chiral centers in α -phellandrene, β -phellandrene, and limonene.
3. Explain how carvone fits the isoprene rule (see essay, "Terpenes and Phenylpropanoids").
4. Using the Cahn–Ingold–Prelog sequence rules, assign priorities to the groups around the chiral carbon in carvone. Draw the structural formulas for (+)- and (–)-carvone with the molecules oriented in the correct position to show the R and S configurations.
5. Explain why limonene elutes from the column before either (+)- or (–)-carvone.
6. Explain why the retention times for both carvone isomers are the same.
7. The toxicity of (+)-carvone in rats is about 400 times greater than that of (–)-carvone. How do you account for this?


The Chemistry of Vision

An interesting and challenging topic for chemists to investigate is how the eye functions. What chemistry is involved in detection of light and transmission of that information to the brain? The first definitive studies on how the eye functions were begun in 1877 by Franz Boll. Boll demonstrated that the red color of the retina of a frog's eye could be bleached yellow by strong light. If the frog was then kept in the dark, the red color of the retina slowly returned. Boll recognized that a bleachable substance had to be connected somehow with the ability of the frog to perceive light.

Most of what is now known about the chemistry of vision is the result of the elegant work of George Wald, Harvard University; his studies, which began in 1933, ultimately resulted in his receiving the Nobel Prize in biology. Wald identified the sequence of chemical events during which light is converted into some form of electrical information that can be transmitted to the brain. Here is a brief outline of that process.

The retina of the eye is made up of two types of photoreceptor cells: **rods** and **cones**. The rods are responsible for vision in dim light, and the cones are responsible for color vision in bright light. The same principles apply to the chemical functioning of the rods and the cones; however, the details of that functioning are less well understood for the cones than for the rods.

Each rod contains several million molecules of **rhodopsin**. Rhodopsin is a complex of a protein, **opsin**, and a molecule derived from Vitamin A, *11-cis*-retinal (sometimes called **retinene**). Little is known about the structure of opsin. The structure of *11-cis*-retinal is shown here.

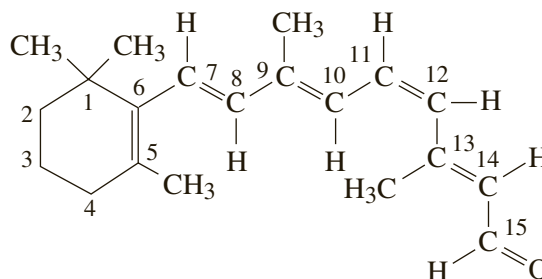


11-*cis*-Retinal

The detection of light involves the initial conversion of *11-cis*-retinal to its all-*trans* isomer. This is the only obvious role of light in this process. The high energy of a quantum of visible light promotes the fission of the π bond between carbons 11 and 12. When the π bond breaks, free rotation about the σ bond in the resulting radical is possible. When the π bond re-forms after such rotation, all-*trans*-retinal results. All-*trans*-retinal is more stable than *11-cis*-retinal, which is the reason the isomerization proceeds spontaneously in the direction shown in the following equation.

The two molecules have different shapes because of their different structures. The 11-*cis*-retinal has a fairly curved shape, and the parts of the molecule on either side of the *cis* double bond tend to lie in different planes. Because proteins have complex and specific three-dimensional shapes (tertiary structures), 11-*cis*-retinal associates with the protein opsin in a particular manner. All-*trans*-retinal has an elongated shape, and the entire molecule tends to lie in a single plane. This different shape for the molecule, compared with that for the 11-*cis* isomer, means that all-*trans*-retinal will have a different association with the protein opsin.

In fact, all-*trans*-retinal associates very weakly with opsin because its shape does not fit the protein. Consequently, the next step after the isomerization of retinal is the dissociation of all-*trans*-retinal from opsin. The opsin protein undergoes a simultaneous change in conformation as the all-*trans*-retinal dissociates.

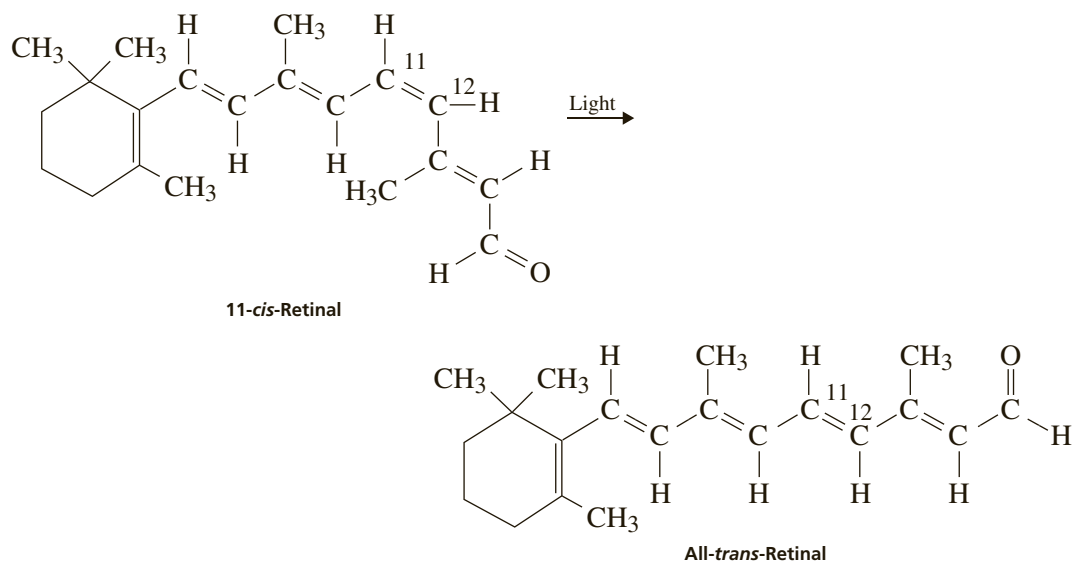

At some time after the 11-*cis*-retinal–opsin complex receives a photon, a message is received by the brain. It was originally thought that either the isomerization of 11-*cis*-retinal to all-*trans*-retinal or the conformational change of the opsin protein was an event that generated the electrical message sent to the brain. Current research,


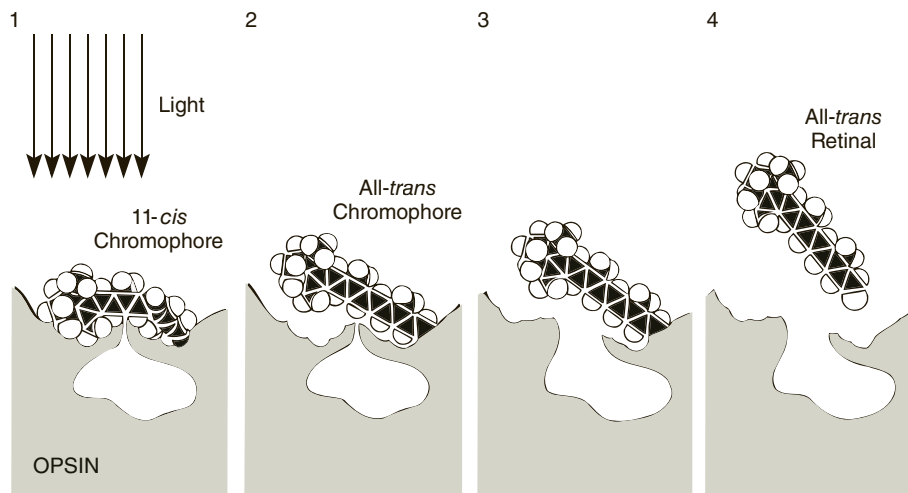



Illustration by Joan Starwood from "Molecular Isomers in Vision," by Ruth Hubbard and Allen Kropf. *Scientific American*, 216:64–76. Reprinted by permission of Joan Starwood.

however, indicates that both these events occur too slowly relative to the speed with which the brain receives the message. Current hypotheses invoke involved quantum mechanical explanations, which regard it as significant that the chromophores (light-absorbing groups) are arranged in a very precise geometrical pattern in the rods and cones, allowing the signal to be transmitted rapidly through space. The main physical and chemical events Wald discovered are illustrated in the figure for easy visualization. The question of how the electrical signal is transmitted still remains unsolved.

Wald was also able to explain the sequence of events by which the rhodopsin molecules are regenerated. After dissociation of all-*trans*-retinal from the protein, the following enzyme-mediated changes occur. All-*trans*-retinal is reduced to the alcohol all-*trans*-retinol, also called all-*trans*-Vitamin A.

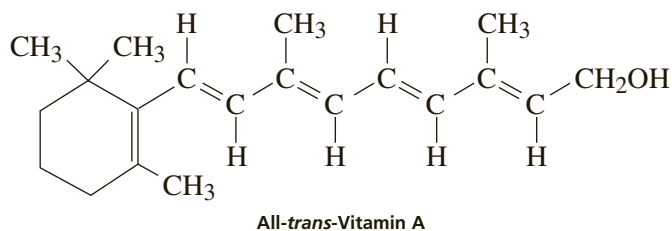


All-*trans*-Vitamin A is then isomerized to its 11-*cis*-Vitamin A isomer. After the isomerization, the 11-*cis*-Vitamin A is oxidized back to 11-*cis*-retinal, which forthwith recombines with the opsin protein to form rhodopsin. The regenerated rhodopsin is then ready to begin the cycle anew, as illustrated in the figure.

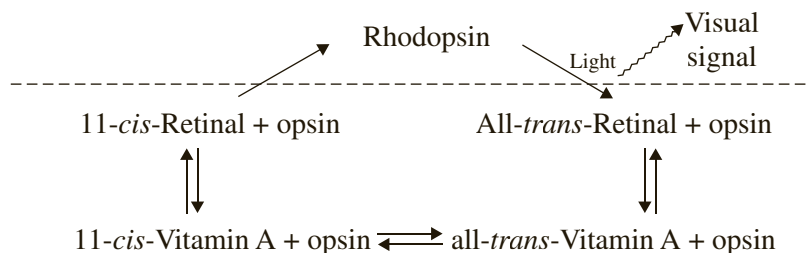

By this process, as little light as 10^{-14} of the number of photons emitted from a typical flashlight bulb can be detected. The conversion of light into isomerized retinal exhibits an extraordinarily high quantum efficiency. Virtually every quantum of light absorbed by a molecule of rhodopsin causes the isomerization of 11-*cis*-retinal to all-*trans*-retinal.

As you can see from the reaction scheme, the retinal derives from Vitamin A, which requires merely the oxidation of a $-\text{CH}_2\text{OH}$ group to a $-\text{CHO}$ group to be converted to retinal. The precursor in the diet that is transformed to Vitamin A is β -carotene. The β -carotene is the yellow pigment of carrots and is an example of a family of long-chain polyenes called **carotenoids**.

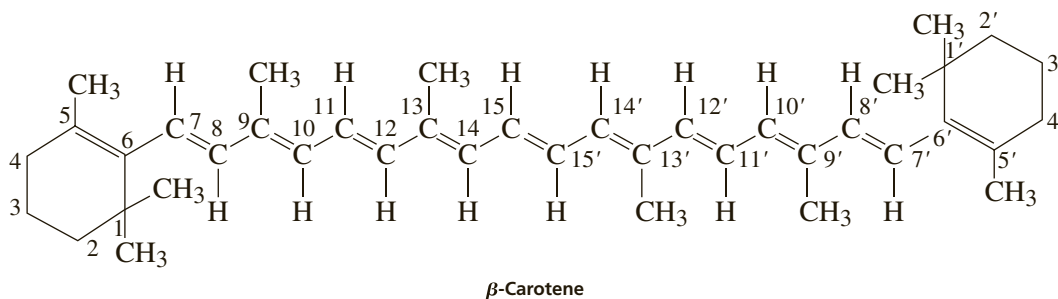


In 1907, Willstätter established the structure of carotene, but it was not known until 1931–1933 that there were actually three isomers of carotene. The α -carotene differs from β -carotene in that the α isomer has a double bond between C_4 and C_5 rather than between C_5 and C_6 , as in the β isomer. The γ isomer has only one ring, identical to the ring in the β isomer, whereas the other ring is opened in the γ form between C_1' and C_6' . The β isomer is by far the most common of the three.

The substance β -carotene is converted to Vitamin A in the liver. Theoretically, one molecule of β -carotene should give rise to two molecules of this vitamin by cleavage of the C_{15} – C_{15}' double bond, but actually only one molecule of Vitamin A is produced from each molecule of carotene. The Vitamin A thus produced is converted to 11-*cis*-retinal within the eye.

Along with the problem of how the electrical signal is transmitted, color perception is also currently under study. In the human eye, there are three kinds of cone cells, which absorb light at 440, 535, and 575 nm, respectively. These cells discriminate among the primary colors. When combinations of them are stimulated, full-color vision is the message received in the brain.

Because all these cone cells use 11-*cis*-retinal as a substrate trigger, it has long been suspected that there must be three different opsin proteins. Recent work has begun to establish how the opsins vary the spectral sensitivity of the cone cells, even though all of them have the same kind of light-absorbing chromophore.

Retinal is an aldehyde, and it binds to the terminal amino group of a lysine residue in the opsin protein to form a Schiff base, or imine linkage ($>C=N-$). This imine linkage is believed to be protonated (with a plus charge) and to be stabilized by being located near a negatively charged amino acid residue of the protein chain. A second negatively charged group is thought to be located near the 11-*cis* double bond. Researchers have recently shown, from synthetic models that use a simpler protein than opsin itself, that forcing these negatively charged groups to be located at different distances from the imine linkage causes the absorption maximum of the 11-*cis*-retinal chromophore to be varied over a wide enough range to explain color vision.

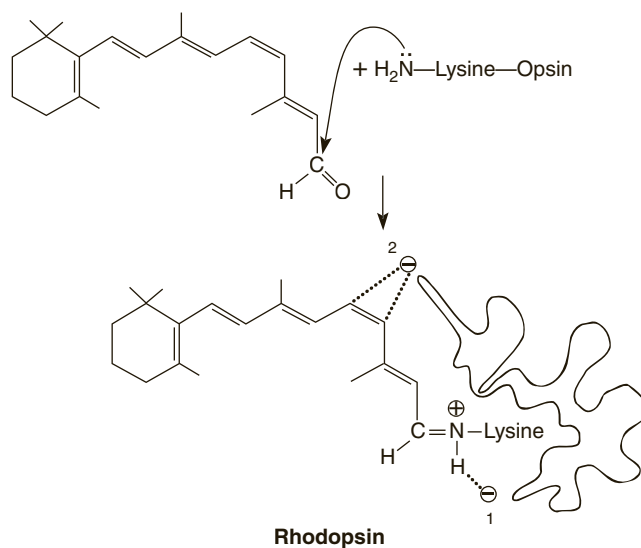

Whether there are actually three different opsin proteins, or whether there are just three different conformations of the same protein in the three types of cone cells, will not be known until further work is completed on the structure of the opsin or opsins.

REFERENCES

- Borman, S. New Light Shed on Mechanism of Human Color Vision. *Chem. Eng. News* **1992**, (Apr 6), 27.
- Fox, J. L. Chemical Model for Color Vision Resolved. *Chem. Eng. News* **1979**, 57 (46) (Nov 12), 25. A review of articles by Honig and Nakanishi in the *J. Am. Chem. Soc.* **1979**, 101, 7082, 7084, 7086.
- Hubbard, R.; and Kropf, A. Molecular Isomers in Vision. *Sci. Am.* **1967**, 216 (Jun), 64.
- Hubbard, R.; and Wald, G. Pauling and Carotenoid Stereochemistry. In *Structural Chemistry and Molecular Biology*; Rich A., and Davidson N., Eds.; W. H. Freeman: San Francisco, 1968.
- MacNichol, E. F., Jr. Three Pigment Color Vision. *Sci. Am.* **1964**, 211 (Dec), 48.
- Model Mechanism May Detail Chemistry of Vision. *Chem. Eng. News* **1985**, 63 (Jan 7), 40.
- Rushton, W. A. H. Visual Pigments in Man. *Sci. Am.* **1962**, 207 (Nov), 120.
- Wald, G. Life and Light. *Sci. Am.* **1959**, 201 (Oct), 92.
- Zurer, P. S. The Chemistry of Vision. *Chem. Eng. News* **1983**, 61 (Nov 28), 24.

Isolation of Chlorophyll and Carotenoid Pigments from Spinach



Isolation of a natural product

Extraction

Column chromatography

Thin-layer chromatography

Photosynthesis in plants takes place in organelles called **chloroplasts**. Chloroplasts contain a number of colored compounds (pigments) that fall into two categories, **chlorophylls** and **carotenoids**.

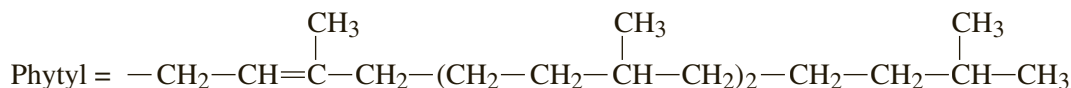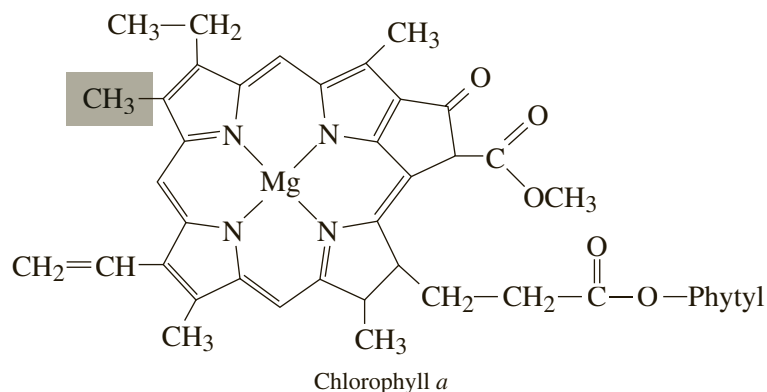

Carotenoids are yellow pigments that are also involved in the photosynthetic process. The structures of α - and β -**carotene** are given in the essay preceding this experiment. In addition, chloroplasts also contain several oxygen-containing derivatives of carotenes, called **xanthophylls**.

In this experiment, you will extract the chlorophyll and carotenoid pigments from spinach leaves using acetone as the solvent. The pigments will be separated by column chromatography using alumina as the adsorbent. Increasingly polar solvents will be used to elute the various components from the column. The colored fractions collected will then be analyzed using thin-layer chromatography. It should be possible for you to identify most of the pigments already discussed on your thin-layer plate after development.

Chlorophylls are the green pigments that act as the principal photoreceptor molecules of plants. They are capable of absorbing certain wavelengths of visible light that are then converted by plants into chemical energy. Two forms of these pigments

found in plants are **chlorophyll a** and **chlorophyll b**. The two forms are identical, except that the methyl group that is shaded in the structural formula of chlorophyll *a* is replaced by a —CHO group in chlorophyll *b*. **Pheophytin a** and **pheophytin b** are identical to chlorophyll *a* and chlorophyll *b*, respectively, except that in each case the magnesium ion Mg^{2+} has been replaced by two hydrogen ions 2H^+ .

REQUIRED READING

Review: Techniques 5 and 6

Technique 7	Reaction Methods, Section 7.9
Technique 12	Extractions, Separations, and Drying Agents, Sections 12.5 and 12.9
Technique 20	Thin-Layer Chromatography

New: Technique 19
 Essay

Column Chromatography
 The Chemistry of Vision

SPECIAL INSTRUCTIONS

Hexane and acetone are both highly flammable. Avoid using flames while working with these solvents. Perform the thin-layer chromatography in the hood. The procedure calls for a centrifuge tube with a tight-fitting cap. If this is not available, you can use a vortex mixer for mixing the liquids. Another alternative is to use a cork to stopper the tube; however, the cork will absorb some liquid.

Fresh spinach is preferable to frozen spinach. Because of handling, frozen spinach contains additional pigments that are difficult to identify. Because the pigments are light-sensitive and can undergo oxidation, you should work quickly. Samples should be stored in closed containers and kept in the dark when possible. The column chromatography procedure takes less than 15 minutes to perform and cannot be stopped until it is completed. It is important, therefore, that you have all the materials needed for this part of the experiment prepared in advance and that you be thoroughly familiar with the procedure before running the column. If you need to prepare the 70% hexane–30% acetone solvent mixture, be sure to mix it thoroughly before using.

SUGGESTED WASTE DISPOSAL

Dispose of all organic solvents in the container for nonhalogenated organic solvents. Place the alumina in the container designated for wet alumina.

NOTES TO THE INSTRUCTOR


The column chromatography should be performed with activated alumina from EM Science (No. AX0612-1). The particle sizes are 80–120 mesh, and the material is Type F-20. Dry the alumina overnight in an oven at 110°C and store it in a tightly sealed bottle. Alumina older than several years may need to be dried for a longer time at a higher temperature. Depending on how dry the alumina is, solvents of different polarity will be required to elute the components from the column.

For thin-layer chromatography, use flexible silica gel plates from Whatman with a fluorescent indicator (No. 4410 222). If the TLC plates have not been purchased recently, place them in an oven at 100°C for 30 minutes and store them in a desiccator until use.

If you use different alumina or different thin-layer plates, try the experiment before using it in class. Materials other than those specified here may give different results than indicated in this experiment.

PROCEDURE

Part A. Extraction of the Pigments

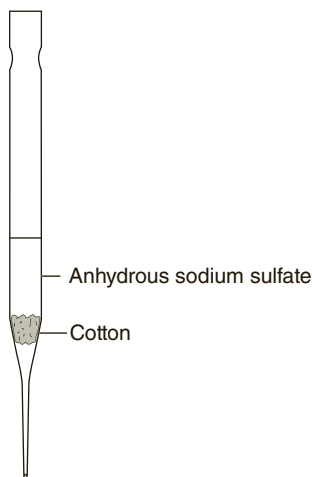

Column for drying extract.

Weigh about 0.5 g of fresh (or 0.25 g of frozen) spinach leaves (avoid using stems or thick veins). Fresh spinach is preferable, if available. If you must use frozen spinach, dry the thawed leaves by pressing them between several layers of paper towels. Cut or tear the spinach leaves into small pieces and place them in a mortar along with 1.0 mL of cold acetone. Grind with a pestle until the spinach leaves have been broken into particles too small to be seen clearly. If too much acetone has evaporated, you may need to add an additional portion of acetone (0.5–1.0 mL) to perform the following step. Using a Pasteur pipette, transfer the mixture to a centrifuge tube. Rinse the mortar and pestle with 1.0 mL of cold acetone and transfer the remaining mixture to the centrifuge tube. Centrifuge the mixture (be sure to balance the centrifuge). Using a Pasteur pipette, transfer the liquid to a centrifuge tube with a tight-fitting cap (see “Special Instructions” if one is not available).

Add 2.0 mL of hexane to the tube, cap the tube, and shake the mixture thoroughly. Then add 2.0 mL of water and shake thoroughly with occasional venting. Centrifuge the mixture to break the emulsion, which usually appears as a cloudy green layer in the middle of the mixture. Remove the bottom aqueous layer with a Pasteur pipette. Using a Pasteur pipette, prepare a column containing anhydrous sodium sulfate to dry the remaining hexane layer, which contains the dissolved pigments. Put a plug of cotton into a Pasteur pipette (5¾-inch) and tamp it into position using a glass rod. The correct position of the cotton is shown in the figure. Add about 0.5 g of powdered or granular anhydrous sodium sulfate and tap the column with your finger to pack the material.

Clamp the column in a vertical position and place a dry test tube (13 × 100-mm) under the bottom of the column. Label this test tube with an *E* for *extract* so that you don't confuse it with the test tubes you will be working with later in this experiment. With a Pasteur pipette, transfer the hexane layer to the column. When all the solution has drained, add 0.5 mL of hexane to the column to extract all the pigments from the drying agent. Evaporate the solvent by placing the test tube in a warm water bath (40–60°C) and directing a stream of nitrogen gas (or dry air) into the test tube. Dissolve the residue in 0.5 mL of hexane. Stopper the test tube and place it in your drawer until you are ready to run the alumina chromatography column.

Part B. Column Chromatography

Introduction

The pigments are separated on a column packed with alumina. Although there are many components in your sample, they usually separate into two main bands on the column. The first band to pass through the column is yellow and consists of the carotenes. This band may be less than 1 mm wide and it may pass through the column rapidly. It is easy to miss seeing the band as it passes through the alumina. The second band consists of all the other pigments discussed in the introduction

to this experiment. Although it consists of both green and yellow pigments, it appears as a green band on the column. The green band spreads out on the column more than the yellow band, and it moves more slowly. Occasionally, the yellow and green components in this band will separate as the band moves down the column. If this begins to occur, you should change to a solvent of higher polarity so that they come out as one band. As the samples elute from the column, collect the yellow band (carotenes) in one test tube and the green band in another test tube.

Because the moisture content of the alumina is difficult to control, different samples of alumina may have different activities. The activity of the alumina is an important factor in determining the polarity of the solvent required to elute each band of pigments. Several solvents with a range of polarities are used in this experiment. The solvents and their relative polarities follow:

Hexane	↓	increasing
70% hexane–30% acetone		polarity
Acetone		
80% acetone–20% methanol	↓	

A solvent of lower polarity elutes the yellow band; a solvent of higher polarity is required to elute the green band. In this procedure, you first try to elute the yellow band with hexane. If the yellow band does not move with hexane, you then add the next more polar solvent. Continue this process until you find a solvent that moves the yellow band. When you find the appropriate solvent, continue using it until the yellow band is eluted from the column. When the yellow band is eluted, change to the next more polar solvent. When you find a solvent that moves the green band, continue using it until the green band is eluted. Remember that occasionally a second yellow band will begin to move down the column before the green band moves. This yellow band will be much wider than the first one. If this occurs, change to a more polar solvent. This should bring all the components in the green band down at the same time.

Advance Preparation

Before running the column, assemble the following glassware and liquids. Obtain five dry test tubes (16 × 100-mm), and number them 1 through 5. Prepare two dry Pasteur pipettes with bulbs attached. Calibrate one of them to deliver a volume of about 0.25 mL. Place 10.0 mL hexane, 6.0 mL 70% hexane–30% acetone solution (by volume), 6.0 mL acetone, and 6.0 mL 80% acetone–20% methanol (by volume) into four separate containers. Clearly label each container.

Prepare a chromatography column packed with alumina. Place a *loose* plug of cotton in a Pasteur pipette (5¾-inch), and push it *gently* into position using a glass rod (see the figure for the correct position of the cotton). Add 1.25 g of alumina (EM Science, No. AX0612-1) to the pipette¹ while tapping the column gently with your finger. When all the alumina has been added, tap the column with your finger for several seconds to ensure that the alumina is tightly packed. Clamp the column in a vertical position so that the bottom of the column is just above the height of the test tubes you will be using to collect the fractions. Place test tube 1 under the column.

¹ As an option, students may prepare a microfunnel from a 1-mL disposable plastic pipette. The microfunnel is prepared by (1) cutting the bulb in half with a scissors and (2) cutting the stem at an angle about ½-inch below the bulb. This funnel can be placed in the top of the column (Pasteur pipette) to aid in filling the column with alumina or with the solvents.

NOTE: Read the following procedure on running the column. The chromatography procedure takes less than 15 minutes, and you cannot stop until all the material is eluted from the column. You must have a good understanding of the whole procedure before running the column.

Running the Column

Using a Pasteur pipette, slowly add about 3.0 mL of hexane to the column. The column must be completely moistened by the solvent. Drain the excess hexane until the level of hexane reaches the top of the alumina. Once you have added hexane to the alumina, the top of the column must not be allowed to run dry. If necessary, add more hexane.

NOTE: It is essential that the liquid level not be allowed to drain below the surface of the alumina at any point during the procedure.

When the level of the hexane reaches the top of the alumina, add about half (0.25 mL) of the dissolved pigments to the column. Leave the remainder in the test tube for the thin-layer chromatography procedure. (Put a stopper on the tube and place it back in your drawer.) Continue collecting the eluent in test tube 1. Just as the pigment solution penetrates the column, add 1 mL of hexane and drain until the surface of the liquid has reached the alumina.

Add about 4 mL of hexane. If the yellow band begins to separate from the green band, continue to add hexane until the yellow band passes through the column. If the yellow band does not separate from the green band, change to the next more polar solvent (70% hexane–30% acetone). When changing solvents, do not add the new solvent until the last solvent has nearly penetrated the alumina. When the appropriate solvent is found, add this solvent until the yellow band passes through the column. Just before the yellow band reaches the bottom of the column, place test tube 2 under the column. When the eluent becomes colorless again (the total volume of the yellow material should be less than 2 mL), place test tube 3 under the column.


Add several mL of the next more polar solvent when the level of the last solvent is almost at the top of the alumina. If the green band moves down the column, continue to add this solvent until the green band is eluted from the column. If the green band does not move or if a diffuse yellow band begins to move, change to the next more polar solvent. Change solvents again if necessary. Collect the green band in test tube 4. When there is little or no green color in the eluent, place test tube 5 under the column and stop the procedure.

Using a warm water bath (40–60°C) and a stream of nitrogen gas, evaporate the solvent from the tube containing the yellow band (tube 2), the tube containing the green band (tube 4), and the tube containing the original pigment solution (tube E). As soon as all the solvent has evaporated from each of the tubes, remove them from the water bath. Do not allow any of the tubes to remain in the water bath after the solvent has evaporated. Stopper the tubes and place them in your drawer.

Part C. Thin-Layer Chromatography

Preparing the TLC Plate

Technique 20 describes the procedures for thin-layer chromatography. Use a 10-cm × 3.3-cm TLC plate (Whatman Silica Gel Plates No. 4410 222). These plates have a flexible backing but should not be bent excessively. Handle them carefully or the adsorbent may flake off them. Also, you should handle them only by the edges; the surface



Preparing the TLC plate.

should not be touched. Using a lead pencil (not a pen), *lightly* draw a line across the plate (short dimension) about 1 cm from the bottom (see figure). Using a centimeter ruler, move its index about 0.6 cm in from the edge of the plate and lightly mark off three 1-cm intervals on the line. These are the points at which the samples will be spotted.

Prepare three micropipettes to spot the plate. The preparation of these pipettes is described and illustrated in Technique 20, Section 20.4. Prepare a TLC development chamber with 70% hexane–30% acetone (see Technique 20, Section 20.5). A beaker covered with aluminum foil or a wide-mouth screw-cap bottle is a suitable container to use (see Technique 20, Figure 20.5). The backing on the TLC plates is thin, so if they touch the filter paper liner of the development chamber *at any point*, solvent will begin to diffuse onto the absorbent surface at that point. To avoid this, be sure that the filter paper liner does not go completely around the inside of the container. A space about 2 inches wide must be provided.

Using a Pasteur pipette, add two drops of 70% hexane–30% acetone to each of the three test tubes containing dried pigments. Swirl the tubes so that the drops of solvent dissolve as much of the pigments as possible. The TLC plate should be spotted with three samples: the extract, the yellow band from the column, and the green band. For each of the three samples, use a different micropipette to spot the sample on the plate. The correct method of spotting a TLC plate is described in Technique 20, Section 20.4. Take up part of the sample in the pipette (don't use a bulb; capillary action will draw up the liquid). For the extract (tube labeled E) and the green band (tube 4), touch the plate once *lightly* and let the solvent evaporate. The spot should be no larger than 2 mm in diameter and should be a fairly dark green. For the yellow band (tube 2), repeat the spotting technique 5–10 times until the spot is a definite yellow. Let the solvent evaporate completely between successive applications and spot the plate in exactly the same position each time. Save the liquid samples in case you need to repeat the TLC.

Developing the TLC Plate

Place the TLC plate in the development chamber, making sure that the plate does not come in contact with the filter paper liner. Remove the plate when the solvent front is 1–2 cm from the top of the plate. Using a lead pencil, mark the position of the solvent front. As soon as the plates have dried, outline the spots with a pencil and indicate the colors. This is important to do soon after the plates have dried because some of the pigments will change color when exposed to the air.

Analysis of the Results

In the crude extract, you should be able to see the following components (in order of decreasing R_f values):

- Carotenes (1 spot) (yellow-orange)
- Pheophytin *a* (gray, may be nearly as intense as chlorophyll *b*)
- Pheophytin *b* (gray, may not be visible)
- Chlorophyll *a* (blue-green, more intense than chlorophyll *b*)
- Chlorophyll *b* (green)
- Xanthophylls (possibly three spots: yellow)

Depending on the spinach sample, the conditions of the experiment, and how much sample was spotted on the TLC plate, you may observe other pigments. These additional components can result from air oxidation, hydrolysis, or other chemical

reactions involving the pigments discussed in this experiment. It is common to observe other pigments in samples of frozen spinach. It is also common to observe components in the green band that were not present in the extract.

Identify as many of the spots in your samples as possible. Determine which pigments were present in the yellow band and which were present in the green band. Draw a picture of the TLC plate in your notebook. Label each spot with its color and its identity, where possible. Calculate the R_f values for each spot produced by chromatography of the extract (see Technique 20, Section 20.9). At the instructor's option, submit the TLC plate with your report.

QUESTIONS

1. Why are the chlorophylls less mobile on column chromatography, and why do they have lower R_f values than the carotenes?
2. Propose structural formulas for pheophytin *a* and pheophytin *b*.
3. What would happen to the R_f values of the pigments if you were to increase the relative concentration of acetone in the developing solvent?
4. Using your results as a guide, comment on the purity of the material in the green and yellow bands; that is, did each band consist of a single component?


Ethanol and Fermentation Chemistry

The fermentation processes involved in making bread, making wine, and brewing are among the oldest chemical arts. Even though fermentation had been known as an art for centuries, it was not until the nineteenth century that chemists began to understand this process from the point of view of science. In 1810, Gay-Lussac discovered the general chemical equation for the breakdown of sugar into ethanol and carbon dioxide. The manner in which the process took place was the subject of much conjecture until Louis Pasteur began his thorough examination of fermentation. Pasteur demonstrated that yeast was required in the fermentation. He was also able to identify other factors that controlled the action of the yeast cells. His results were published in 1857 and 1866.


For many years, scientists believed that the transformation of sugar into ethanol and carbon dioxide by yeasts was inseparably connected with the life process of the yeast cell. This view was abandoned in 1897, when Büchner demonstrated that yeast extract would bring about alcoholic fermentation in the absence of any yeast cells. The fermenting activity of yeast is due to a remarkably active catalyst of biochemical origin, the enzyme zymase. It is now recognized that most of the chemical transformations that occur in living cells of plants and animals are brought about by enzymes. These enzymes are organic compounds, generally proteins, and establishment of the structures and reaction mechanisms of these compounds is an active field of present-day research. Zymase is now known to be a complex of at least 22 separate enzymes, each of which catalyzes a specific step in the fermentation reaction sequence.

Enzymes display an extraordinary specificity—a given enzyme acts on a specific compound or a closely related group of compounds. Thus, zymase acts on only a few select sugars and not on all carbohydrates; the digestive enzymes of the alimentary tract are equally specific in their activity.

The chief sources of sugars for fermentation are the various starches and the molasses residue obtained from refining sugar. Corn (maize) is the chief source of starch in the United States, and ethyl alcohol made from corn is commonly known as **grain alcohol**. In preparing alcohol from corn, the grain, with or without the germ, is ground and cooked to give the **mash**. The enzyme diastase is added in the form of **malt** (sprouted barley that has been dried in air and ground to a powder) or of a mold such as *Aspergillus oryzae*. The mixture is kept at 40°C until all the starch has been converted to the sugar **maltose** by hydrolysis of ether and acetal bonds. This solution is known as the **wort**.

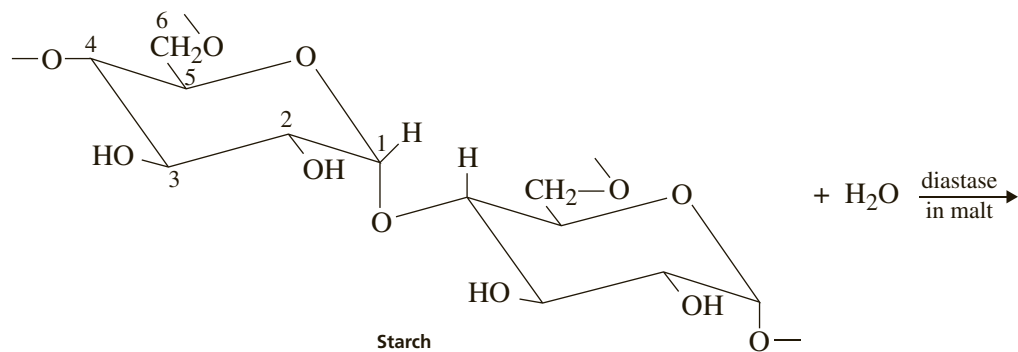

This is a glucose polymer with 1,4- and 1,6- glycosidic linkages. The linkages at C-1 are α .

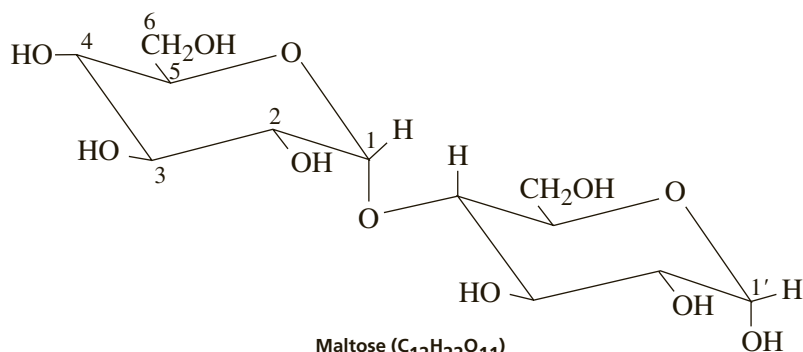


The α linkage still exists at C-1.
The —OH is shown α at the 1' position (axial), but it can also be β (equatorial).

The wort is cooled to 20°C and diluted with water to 10% maltose, and a pure yeast culture is added. The yeast culture is usually a strain of *Saccharomyces cerevisiae* (or *ellipsoidus*). The yeast cells secrete two enzyme systems: maltase, which converts the maltose into glucose, and zymase, which converts the glucose into carbon dioxide and alcohol. Heat is liberated, and the temperature must be kept below 35°C by cooling to prevent destruction of the enzymes. Oxygen in large amounts is initially necessary for the optimum reproduction of yeast cells, but the actual production of alcohol is anaerobic. During fermentation, the evolution of carbon dioxide soon establishes anaerobic conditions. If oxygen were freely available, only carbon dioxide and water would be produced.

After 40–60 hours, fermentation is complete, and the product is distilled to remove the alcohol from solid matter. The distillate is fractionated by means of an efficient column. A small amount of acetaldehyde (bp 21°C) distills first and is followed by 95% alcohol. Fusel oil is contained in the higher-boiling fractions. The fusel oil consists of a mixture of higher alcohols, chiefly 1-propanol, 2-methyl-1-propanol, 3-methyl-1-butanol, and 2-methyl-1-butanol. The exact composition of fusel oil varies considerably; it particularly depends on the type of raw material that is fermented. These higher alcohols are not formed by fermentation of glucose. They arise from certain amino acids derived from the proteins present in the raw material and the yeast. These fusel oils cause the headaches associated with drinking alcoholic beverages.

Industrial alcohol is ethyl alcohol used for nonbeverage purposes. Most commercial alcohol is denatured to avoid payment of taxes, the biggest cost in the price of liquor. The denaturants render the alcohol unfit for drinking. Methanol, aviation fuel, and

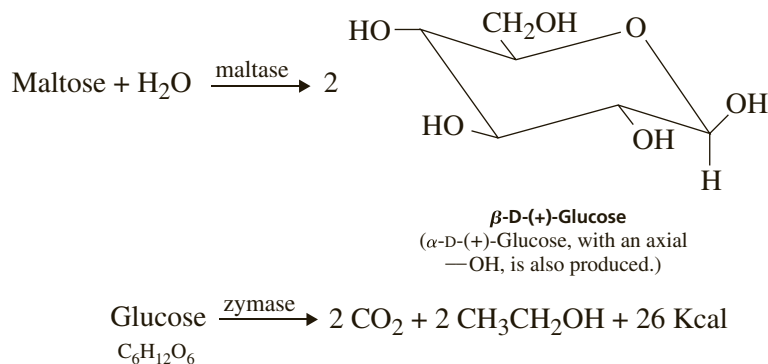

other substances are used for this purpose. The difference in price between taxed and nontaxed alcohol is more than \$20 a gallon. Before efficient synthetic processes were developed, the chief source of industrial alcohol was fermented blackstrap molasses, the noncrystallizable residue from refining cane sugar (sucrose). Most industrial ethanol in the United States is now manufactured from ethylene, a product of the “cracking” of petroleum hydrocarbons. By reaction with concentrated sulfuric acid, ethylene becomes ethyl hydrogen sulfate, which is hydrolyzed to ethanol by dilution with water. The alcohols 2-propanol, 2-butanol, 2-methyl-2-propanol, and higher secondary and tertiary alcohols are also produced on a large scale from alkenes derived from cracking.

Yeasts, molds, and bacteria are used commercially for the large-scale production of various organic compounds. An important example, in addition to ethanol production, is the anaerobic fermentation of starch by certain bacteria to yield 1-butanol, acetone, ethanol, carbon dioxide, and hydrogen.

For additional information on the production of ethanol, see the essay Biofuels that precedes Experiment 28. In this essay, the production of ethanol from corn for use in automobiles is discussed, along with the production of ethanol from other sources such as plant cellulose.

REFERENCES

- Amerine, M. A. Wine. *Sci. Am.* **1964**, 211 (Aug), 46.
- Hallberg, D. E. Fermentation Ethanol. *ChemTech* **1984**, 14 (May), 308.
- Ough, C. S. Chemicals Used in Making Wine. *Chem. Eng. News* **1987**, 65 (Jan 5), 19.
- Van Koevering, T. E.; Morgan, M. D.; and Younk, T. J. The Energy Relationships of Corn Production and Alcohol Fermentation. *J. Chem. Educ.* **1987**, 64 (Jan), 11.
- Webb, A. D. The Science of Making Wine. *Am. Sci.* **1984**, 72 (Jul–Aug), 360.
- Students wanting to investigate alcoholism and possible chemical explanations for alcohol addiction may consult the following references (note that these ideas are mainly hypotheses):
- Cohen, G.; and Collins, M. Alkaloids from Catecholamines in Adrenal Tissue: Possible Role in Alcoholism. *Science* **1970**, 167, 1749.
- Davis, V. E.; and Walsh, M. J. Alcohol Addiction and Tetrahydropapaveroline. *Science* **1970**, 169, 1105.
- Davis, V. E.; and Walsh, M. J. Alcohols, Amines, and Alkaloids: A Possible Biochemical Basis for Alcohol Addiction. *Science* **1970**, 167, 1005.
- Seevers, M. H.; Davis, V. E.; and Walsh, M. J. Morphine and Ethanol Physical Dependence: A Critique of a Hypothesis. *Science* **1970**, 170, 1113.
- Yamanaka, Y.; Walsh, M. J.; and Davis, V. E. Salsolinol, an Alkaloid Derivative of Dopamine Formed in Vitro during Alcohol Metabolism. *Nature* **1970**, 227, 1143.


Ethanol from Sucrose

Fermentation

Fractional distillation

Azeotropes

Either sucrose or maltose can be used as the starting material for making ethanol. Sucrose is a disaccharide with the formula $C_{12}H_{22}O_{11}$. It has one glucose molecule combined with fructose. Maltose consists of two glucose molecules. The enzyme **invertase** is used to catalyze the hydrolysis of sucrose. **Maltase** is more effective in catalyzing the hydrolysis of maltose. The hydrolysis of maltose is discussed in the essay on ethanol and fermentation. **Zymase** is used to convert the hydrolyzed sugars to alcohol and carbon dioxide. Pasteur observed that growth and fermentation were promoted by adding small amounts of mineral salts to the nutrient medium. Later, it was found that before fermentation actually begins, the hexose sugars combine with phosphoric acid, and the resulting hexose-phosphoric acid combination is then degraded into carbon dioxide and ethanol. The carbon dioxide is not wasted in the commercial process, because it is converted to dry ice.

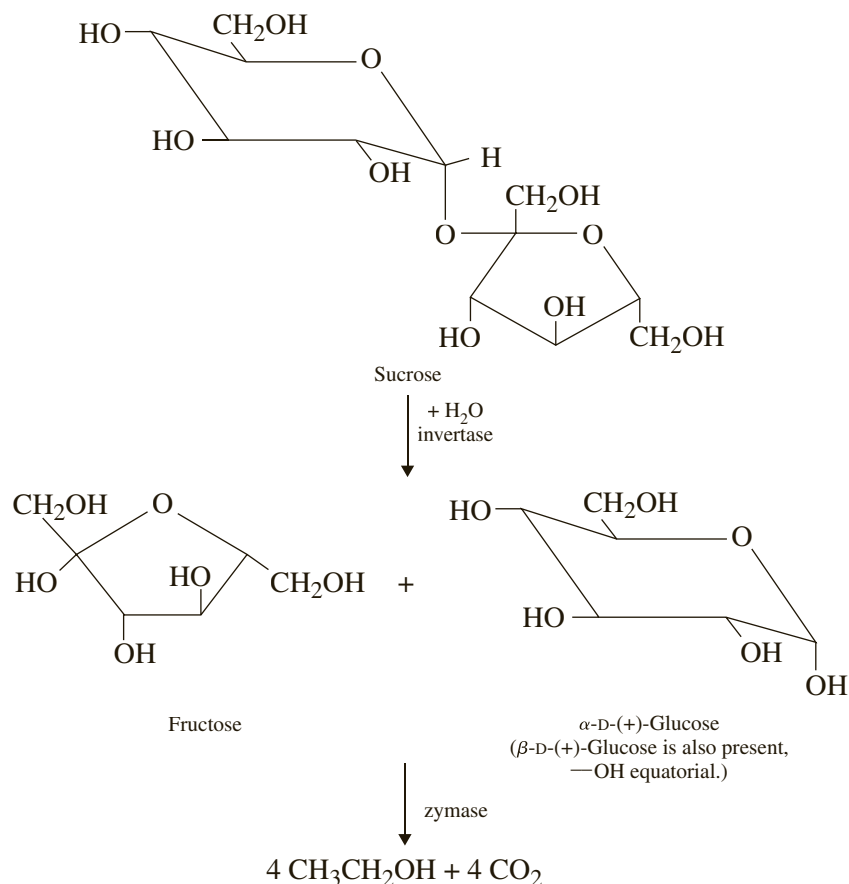

The fermentation is inhibited by its end product ethanol; it is not possible to prepare solutions containing more than 10–15% ethanol by this method. More concentrated ethanol can be isolated by fractional distillation. Ethanol and water form an azeotropic mixture consisting of 95% ethanol and 5% water by weight, which is the most concentrated ethanol that can be obtained by fractionation of dilute ethanol–water mixtures.

REQUIRED READING

Review: Technique 8 Sections 8.3 and 8.4
 Technique 13 Physical Constants of Liquids

New: Technique 15 Fractional Distillation, Azeotropes
 Essay Ethanol and Fermentation Chemistry

SPECIAL INSTRUCTIONS

Start the fermentation at least 1 week before the period in which the ethanol will be isolated. When the aqueous ethanol solution is to be separated from the yeast cells, it is important to transfer carefully as much of the clear, supernatant liquid as possible without agitating the mixture.

NOTES TO THE INSTRUCTOR

Because the volume of the fermentation mixture is only about 20 mL, it is necessary to use an external heat source to maintain a temperature of 30–35°C. An incubator will provide the necessary temperature control. One can make a simple incubator by placing a cardboard box over a light bulb that is turned on during the fermentation. Be sure that the box does not touch the light bulb and has adequate clearance. Use aluminum foil to seal any openings and to help reflect the heat inward.

Using Saran™ or other plastic wrap to cover the opening of the flask is usually sufficient to exclude oxygen from the fermentation mixture.

One method of insulating the air condenser used for the fractional distillation column is provided by using two layers of clear flexible tubing (PVC) over the air condenser. For a ½-inch diameter column, use ½-inch I.D. × ⅝-inch O.D. plastic tubing on the inside and ⅝-inch I.D. × ⅞-inch O.D. tubing on the outside. Cut the tubing into 3½-inch lengths. Make a slit from end to end so that they can slip over the column. Slit the tubing using a sharp scissors or a razor knife with a proper handle.

CAUTION



Do not use a razor blade or you may get badly cut.

The clear tubing allows you to see what is going on in the column and also provides some insulation. Another method of insulating the fractionating column is to

wrap the air condenser with a cotton pad about 3 inches square. Prepare the cotton pad by covering both sides of one layer of cotton with aluminum foil. Wrap this entirely with duct tape to hold the cotton in place and to make a more durable pad.

A convenient, safe, and accurate way to monitor the temperature during the distillation is to use a Vernier LabQuest device with a stainless steel probe or a Vernier LabPro interface with a laptop computer and stainless steel probe (see Technique 13, Section 13.4, and Technique 14, Section 14.5, and Figure 14.12). With both of these methods, students observe a graph of time vs. temperature. Another convenient method is to use a digital thermometer with a stainless steel probe (see Technique 14, Section 14.5, and Figure 14.12); however, most of these devices do not provide a graph of the temperature. All of these methods are more accurate than using non-mercury thermometers. If a glass thermometer is used, the temperature will be most accurate if a partial immersion mercury thermometer is used. See the Instructor's Manual for additional comments about the use of these devices, including suitable stainless steel probes for this experiment.

This experiment can also be performed without doing the fermentation. Provide each student with 20 mL of a 10% ethanol solution. This solution is used in place of the fermentation mixture in the Fractional Distillation section of the Procedure.

PROCEDURE

Fermentation

Place 2.00 g of sucrose in a 50-mL Erlenmeyer flask. Add 18.0 mL of water warmed to 25–30°C, 2.0 mL of Pasteur's salts,¹ and 0.2 g of *dried* baker's yeast. Shake the contents vigorously to mix them, and then cover the opening of the flask with Saran™ wrap tightly secured in place with a rubber band. Oxygen from the atmosphere is excluded from the chemical reaction by this technique. If oxygen were allowed to continue in contact with the fermenting solution, the ethanol could be further oxidized to acetic acid or even all the way to carbon dioxide and water. As long as carbon dioxide continues to be liberated, ethanol is being formed.

Allow the mixture to stand at about 30–35°C for 1 week. After this time, *carefully* move the flask away from the heat source and remove the balloon. Without disturbing the sediment, transfer the clear, supernatant liquid solution to another container with a Pasteur pipette. Try to avoid drawing any of the sediment into the pipette.

If it is not possible to remove the solution completely without drawing up sediment, remove the sediment by centrifugation. Pour equal amounts of the liquid into two centrifuge tubes. After centrifugation for several minutes, decant the liquid away from the solid into another container. The liquid contains ethanol in water, plus smaller amounts of dissolved metabolites (fusel oils) from the yeast. The mixture will be subjected to fractional distillation.

¹ A solution of Pasteur's salts consists of potassium dihydrogen phosphate, 1.0 g; calcium phosphate (monobasic), 0.10 g; magnesium sulfate, 0.10 g; and ammonium tartrate (diammonium salt), 5.0 g, dissolved in 430 mL water.

Fractional Distillation

Assemble the apparatus shown in Technique 15, Figure 15.2; use a 20-mL or 25-mL round-bottom flask in place of the 10-mL flask. Use an aluminum block, if available, for the heat source. Pack the air condenser *uniformly* with about 1 g of stainless steel cleaning pad material. (Use a pad that does not contain soap.)

CAUTION



You should wear heavy cotton gloves when handling the stainless steel cleaning pad. The edges are sharp and can easily cut into the skin.

Wrap the glass section of the air condenser between the two plastic caps with plastic tubing as described in Notes to the Instructor. Alternatively, use the method with a cotton pad (see Notes to the Instructor). Hold the pad in place with tape or twist ties. Place a boiling stone and the fermentation mixture in the round-bottom flask. If a sand bath is used, the apparatus should be clamped so that the bottom half of the flask is buried in the sand.

The temperature during the distillation may be monitored either with a thermometer or a stainless steel temperature probe. If a stainless steel probe is used, it must be used in conjunction with either a digital thermometer or one of the Vernier devices (see Technique 13, Section 13.4 and Technique 14, Figure 14.12). Your instructor will provide instructions about the method that you will use. Insert the thermometer (or probe) so that the bulb is level with or slightly below the cap connecting the Hickman head to the air condenser. Also use a thermometer to monitor the temperature of the heat source. Cover the top of the sand bath (if used) with a square of aluminum foil with a tear from the center of one edge to the middle.

It is important to distill the liquid *slowly* through the fractionating column to get the best possible separation. This can be done by carefully following these instructions: Adjust the temperature of the heat source to achieve a rapid boiling rate in the flask. This will initially require a high setting on the hot plate. When the liquid begins boiling, it is best to turn the heat down immediately and then gradually raise it so that the heat setting required to maintain boiling is at the lowest possible setting. It may be necessary to increase the temperature of the heat source as the distillation proceeds, especially if boiling in the flask stops or if the distillation seems to be taking a very long time. On the other hand, if a lot of liquid quickly fills the column, remove the heat source for a short time so that the liquid drains back into the flask. Once ethanol reaches the top of the column, the temperature in the distillation head will increase to about 78°C and remain at this temperature until the ethanol fraction is distilled.

As distillate condenses in the Hickman head, transfer the liquid from the reservoir to a preweighed 3-mL conical vial. If your Hickman head does not have a side port, it will be necessary to use a 9-inch Pasteur pipette. In the latter case, it is helpful to bend the tip of the pipette slightly by heating it in a flame. The distillate can then be removed without removing the thermometer. Be sure to cap the conical vial used for storage each time after you transfer the distillate. Continue to distill the mixture, and transfer the distillate to the vial until the temperature in the Hickman head increases above 78°C or until the temperature in the Hickman head drops several degrees below 78°C and remains at this lower temperature for 10 minutes or more. You should collect about 0.4 mL of distillate. The distillation should then be interrupted by removing the apparatus from the heat source.

Analysis of Distillate

Determine the total weight of the distillate. Determine the approximate density of the distillate by transferring a known volume of the liquid with an automatic pipette or graduated pipette to a tared vial. Reweigh the vial and calculate the density. This method is good to two significant figures. Using the following table, determine the percentage composition by weight of ethanol in your distillate from the density of your sample. The extent of purification of the ethanol is limited because ethanol and water form a constant-boiling mixture, an azeotrope, with a composition of 95% ethanol and 5% water.

Percentage Ethanol by Weight	Density at 20°C (g/mL)	Density at 25°C (g/mL)
75	0.856	0.851
80	0.843	0.839
85	0.831	0.827
90	0.818	0.814
95	0.804	0.800
100	0.789	0.785

Calculate the percentage yield of alcohol. At the option of the instructor, determine the boiling point of the distillate using a microboiling-point method (Technique 13, Section 13.2). The boiling point of the azeotrope is 78.1°C. Submit the ethanol to the instructor in a labeled vial.²

QUESTIONS

1. Write a balanced equation for the conversion of sucrose into ethanol.
2. By doing some library research, see whether you can find the commercial method or methods used to produce absolute ethanol.
3. Why is the balloon necessary in the fermentation?
4. How does acetaldehyde impurity arise in the fermentation?
5. Diethylacetal can be detected by gas chromatography. How does this impurity arise in fermentation?
6. Calculate how many milliliters of carbon dioxide would be produced theoretically from 2.0 g of sucrose at 25°C and 1 atmosphere pressure.

² A careful analysis by flame-ionization gas chromatography on a typical student prepared ethanol sample provided the following results:

Acetaldehyde	0.060%
Diethylacetal of acetaldehyde	0.005%
Ethanol	88.3% (by hydrometer)
1-Propanol	0.031%
2-Methyl-1-propanol	0.092%
5-Carbon and higher alcohols	0.140%
Methanol	0.040%
Water	11.3% (by difference)

Introduction to Molecular Modeling

Molecular Modeling and Molecular Mechanics

Since the beginnings of organic chemistry, somewhere in the middle of the nineteenth century, chemists have sought to visualize the three-dimensional characteristics of the all-but-invisible molecules that participate in chemical reactions. Concrete models that could be held in the hand were developed. Many kinds of model sets, such as framework, ball-and-stick, and space-filling models, were devised to allow people to visualize the spatial and directional relationships within molecules. These handheld models were interactive, and they could be readily manipulated in space.

Today we can also use the computer to help us visualize these molecules. The computer images are also completely interactive, allowing us to rotate, scale, and change the type of model viewed at the press of a button or the click of a mouse. In addition, the computer can rapidly calculate many properties of the molecules that we view. This combination of visualization and calculation is often called **computational chemistry** or, more colloquially, **molecular modeling**.

Two distinct methods of molecular modeling are commonly used by organic chemists today. The first of these is **quantum mechanics**, which involves the calculation of orbitals and their energies using solutions of the Schrödinger equation. The second method is not based on orbitals at all but instead is founded on our knowledge of the way in which the bonds and angles in a molecule behave. Classical equations that describe the stretching of bonds and the bending of angles are used. This second approach is called **molecular mechanics**. The two types of calculation are used for different purposes and do not calculate the same types of molecular properties. In this essay, molecular mechanics will be discussed.

MOLECULAR MECHANICS

Molecular mechanics (MM) was first developed in the early 1970s by two groups of chemical researchers: the Engler, Andose, and Schleyer group and the Allinger group. In molecular mechanics, a mechanical **force field** is defined that is used to calculate an energy for the molecule under study. The energy calculated is often called the **strain energy** or **steric energy** of the molecule. The force field comprises several components, such as bond-stretching energy, angle-bending energy, and bond-torsion energy. A typical force field expression might be represented by the following composite expression:¹

$$E_{\text{strain}} = E_{\text{stretch}} + E_{\text{angle}} + E_{\text{torsion}} + E_{\text{oop}} + E_{\text{vdW}} + E_{\text{dipole}}$$







¹ Other force fields may be found that include more terms than this one and that contain more sophisticated calculational methods than those shown here.

To calculate the final strain energy for a molecule, the computer systematically changes every bond length, bond angle, and torsional angle in the molecule, recalculating the strain energy each time, keeping each change that minimizes the total energy, and rejecting those that increase the energy. In other words, all the bond lengths and angles are changed until the energy of the molecule is *minimized*.

Each term contained in the composite expression (E_{strain}) is defined in Table 1. All of these terms come from classical physics, not quantum mechanics. We will not discuss every term, but will take E_{stretch} as an illustrative example.

Classical mechanics says that a bond behaves like a spring. Each type of bond in a molecule can be assigned a normal bond length, x_0 . If the bond is stretched or compressed, its potential energy will increase, and there will be a restoring force

Table 1 Some of the Factors Contributing to a Molecular Force Field

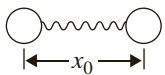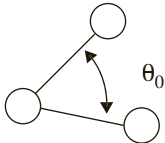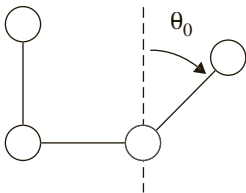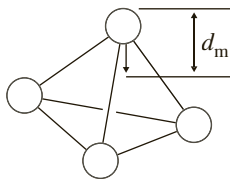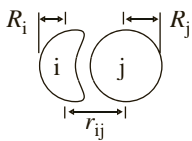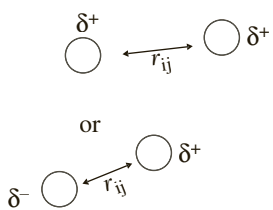Type of Contribution	Illustration	Typical Equation
E_{stretch} (bond stretching)		$E_{\text{stretch}} = \sum_{i=1}^{n_{\text{bonds}}} (k_i/2)(x_i - x_0)^2$
E_{angle} (angle bending)		$E_{\text{angle}} = \sum_{j=1}^{n_{\text{angles}}} (k_j/2)(\theta_j - \theta_0)^2$
E_{torsion} (bond torsion)		$E_{\text{torsion}} = \sum_{k=1}^{n_{\text{torsions}}} (k_k/2)[1 + \text{sp}_k(\cos p_k \theta)]$
E_{oop} (out of plane bending)		$E_{\text{oop}} = \sum_{m=1}^{n_{\text{oops}}} (k_m/2)d_m^2$
E_{vdW} (van der Waals repulsion)		$E_{\text{vdW}} = \sum_{i=1}^{n_{\text{atoms}}} \sum_{j=1}^{n_{\text{atoms}}} (E_i E_j)^{1/2} \left[\frac{1}{a_{ij}^{12}} - \frac{2}{a_{ij}^6} \right]$ $a_{ij} = r_{ij}/(R_i + R_j)$
E_{dipole} (electric dipole repulsion or attraction)		$E_{\text{vdW}} = K \sum_{i=1}^{n_{\text{atoms}}} \sum_{j=i+1}^{n_{\text{atoms}}} Q_i Q_j / r_{ij}^2$

Note: The factors selected here are similar to those in the “Tripos force field” used in the Alchemy III molecular modeling program.

that attempts to restore the bond to its normal length. According to Hooke's Law, the restoring force is proportional to the size of the displacement:

$$F = -k_i (x_1 - x_0) \text{ or } F = -k_i(\Delta x)$$

where k_i is the **force constant** of the bond being studied (that is, the "stiffness" of the spring) and Δx is the change in bond length from the bond's normal length (x_0). The actual energy term that is minimized is given in Table 1. This equation indicates that all the bonds in the molecule contribute to the strain; it is a sum (Σ) starting with the first bond's contribution ($n = 1$) and proceeding through the contributions of all the other bonds (n_{bonds}).


These calculations are based on empirical data. To perform these calculations, the system must be **parameterized** with experimental data. To parameterize, a table of the normal bond lengths (x_0) and force constants (k_i) for every type of bond in the molecule must be created. The program uses these experimental parameters to perform its calculations. The quality of the results from any molecular mechanics approach directly depends on how well the parameterization has been performed for each type of atom and bond that has to be considered. The MM procedure requires each of the factors in Table 1 to have its own parameter table.

Each of the first four terms in Table 1 is treated as a spring in the same manner as discussed for bond stretching. For instance, an angle also has a force constant k that resists a change in the size of the angle θ . In effect, in the first four terms the molecule is treated as a collection of interacting springs, and the energy of this collection of springs must be minimized. In contrast, the last two terms are based on electrostatic or "coulomb" repulsions. Without describing these terms in detail, it should be understood that they must also be minimized.

MINIMIZATION AND CONFORMATION

The object of minimizing the strain energy is to find the lowest energy *conformation* of a molecule. Molecular mechanics does a very good job of finding conformations, because it varies bond distances, bond angles, torsional angles, and the positions of atoms in space. However, most minimizers have some limitations of which users must be aware. Many of the programs use a minimization procedure that will locate a local minimum in the energy but will not necessarily find a global minimum. The figure "Global and local energy minima" that is shown below illustrates the problem.

In the figure, the molecule under consideration has two conformations that represent energy minima for the molecule. Many minimizers will not automatically find the lowest energy conformation, the **global minimum**. The global minimum will be found only when the structure of the starting molecule is already close to the global minimum's conformation. For instance, if the starting structure corresponds to point B on the curve in the figure, then the global minimum will be found. However, if the starting molecule is not close to the global minimum in structure, a **local minimum** (one nearby) may be found. In the figure, if the starting structure corresponds to point A, then a local minimum will be found, instead of the global minimum. Some of the more expensive programs always find the global minimum because they use more sophisticated minimization procedures that depend on random (Monte Carlo) changes instead of sequential ones. However, unless the program has specifically dealt with this problem, the user must be careful to avoid finding a false local minimum when the global minimum is expected. It may be necessary to use several different starting structures to discover the global minimum for a given molecule.

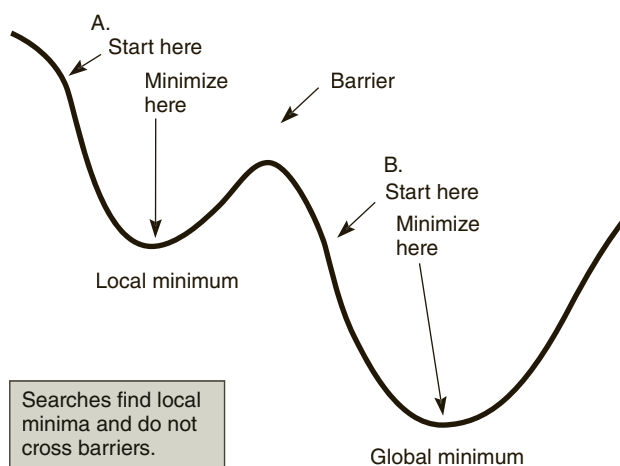

Global and local energy minima.

LIMITATIONS OF MOLECULAR MECHANICS

From our discussion thus far, it should be obvious that molecular mechanics was developed to find the lowest energy conformation of a given molecule or to compare the energies of several conformations of the same molecule. Molecular mechanics calculates a “strain energy,” not a thermodynamic energy such as a heat of formation. Procedures based on quantum mechanics and statistical mechanics are required to calculate thermodynamic energies. Therefore, it is very dangerous to compare the strain energies of two *different* molecules. For instance, molecular mechanics can make a good evaluation of the relative energies of *anti*- and *gauche*-butane conformations, but it cannot fruitfully compare butane and cyclobutane. Isomers can be compared only if they are very closely related. The *cis*- and *trans*-isomers of 1,2-dimethylcyclohexane, or those of 2-butene, can be compared. However, the isomers 1-butene and 2-butene cannot be compared; one is a mono-substituted alkene, whereas the other is disubstituted.

Molecular mechanics will perform the following tasks quite well:

1. It will give good estimates for the actual bond lengths and angles in a molecule.
2. It will find the best conformation for a molecule, but you must watch out for local minima!

Molecular mechanics will not calculate the following properties:

1. It will not calculate thermodynamic properties such as the heat of formation of a molecule.²
2. It will not calculate electron distributions, charges, or dipole moments.
3. It will not calculate molecular orbitals or their energies.
4. It will not calculate infrared, NMR, or ultraviolet spectra.

CURRENT IMPLEMENTATIONS

With time, the most popular version of molecular mechanics has become that developed by Norman Allinger and his research group. The original program from this group was called MM1. The program has undergone constant revisions and improvements, and the current Allinger versions are now designated MM2 and MM3.

²Some of the latest versions are now parameterized to give heats of formation.

However, many other versions of molecular mechanics are now available from both private and commercial sources. Some popular commercial programs that now incorporate their own force fields and parameters include HyperChem, Student Hyperchem, PC Model, MacroModel, Spartan '14, Spartan Student, and Sybyl. You should also realize, however, that many modeling programs do not have molecular mechanics or minimization. These programs will “clean up” a structure that you create by attempting to make every bond length and angle “perfect.” With these programs, every sp^3 carbon will have 109° angles, and every sp^2 carbon will have perfect 120° angles. Using one of these programs is equivalent to using a standard model set that has connectors and bonds with perfect angles and lengths. If you intend to find a molecule's preferred conformation, be sure you use a program that has a force field and performs a true minimization procedure. Also remember that you may have to control the starting structure's geometry in order to find the correct result.

REFERENCES

- Casanova, J. Computer-Based Molecular Modeling in the Curriculum. Computer Series 155. *J. Chem. Educ.* **1993**, 70 (Nov), 904.
- Clark, T. *A Handbook of Computational Chemistry—A Practical Guide to Chemical Structure and Energy Calculations*; John Wiley & Sons: New York, 1985.
- Hypercube, Inc. (www.hyper.com), Gainesville, FL: Hyperchem, Student Hyperchem, and Pocket Hyperchem.
- Lipkowitz, K. B. Molecular Modeling in Organic Chemistry—Correlating Odors with Molecular Structure. *J. Chem. Educ.* **1989**, 66 (Apr), 275.
- Tripos Associates. *Alchemy III—User's Guide*; Tripos Associates: St. Louis, 1992.
- Ulrich, B.; Allinger, N. L. *Molecular Mechanics, ACS Monograph 177*; American Chemical Society: Washington, DC, 1982.
- Wavefunction, Inc. (www.wavefun.com), Irvine, CA: Spartan '14, Spartan Student, and iSpartan. These programs will perform molecular mechanics. iSpartan is a server-connected app for iPod devices.

An Introduction to Molecular Modeling

Molecular modeling
Molecular mechanics

REQUIRED READING

Review: The sections of your lecture textbook dealing with

1. Conformation of cyclic and acyclic compounds
2. The energies of alkenes with respect to degree of substitution
3. The relative energies of *cis*- and *trans*-alkenes

New: Essay Molecular Modeling and Molecular Mechanics

SPECIAL INSTRUCTIONS

To perform this experiment, you must use computer software that has the ability to perform molecular mechanics (MM2 or MM3) calculations with minimization of the strain energy. Either your instructor will provide directions for using the software or you will be given a handout with instructions.

NOTES TO THE INSTRUCTOR

This molecular mechanics experiment was devised using the modeling program Spartan; however, it should be possible to use many other implementations of molecular mechanics. Some of the other capable programs available are HyperChem, Student HyperChem, PCModel, Macromodel, Spartan '14, Spartan Student, and Sybyl. You will have to provide your students with an introduction to your specific implementation. The introduction should show students how to build a molecule, how to minimize its energy, and how to load and save files. Students will also need to be able to measure bond lengths and bond angles.

19A EXPERIMENT 19A

The Conformations of n-Butane: Local Minima

The acyclic butane molecule has several conformations derived by rotation about the C2–C3 bond. The relative energies of these conformations have been well established experimentally and are listed in the following table.

Conformation	Torsional Angle	Relative Energy (kcal/mol)	Relative Energy (kJ/mol)	Types of Strain
Syn	0°	6.0	25.0	Steric/ torsional
Gauche	60°	1.0	4.2	Steric
Eclipsed-120	120°	3.4	14.2	Torsional
Anti	180°	0	0	No strain

In this section, we will show that although molecular mechanics does not calculate the precise thermodynamic energies for the conformations of butane, it will give strain energies that predict the *order* of stability correctly. We will also investigate the difference between a local minimum and a global minimum.

When you construct a butane skeleton, you might expect the minimizer to always arrive at the *anti* conformation (lowest energy). In fact, for most molecular mechanics programs, this will happen only if you bias the minimizer by starting with a butane skeleton that closely resembles the *anti* conformation. If this is done, the minimizer will find the *anti* conformation (the *global* minimum). However, if a skeleton is constructed that does not closely resemble the *anti* conformation, the butane will usually minimize to the *gauche* conformation (the nearest *local* minimum) and not proceed to the global minimum. For the two staggered conformations, you will begin by constructing your starting butane molecules with torsional angles slightly removed from the two minima. The eclipsed conformations, however, will be set on the exact angles to see if they will minimize. Your data should be recorded in a table with the following headings: *Starting Angle*, *Minimized Angle*, *Final Conformation*, and *Minimized Energy*.

Your program should have a feature that allows you to set bond lengths, bond angles, and torsional angles.¹ If it does, you can merely select the torsion angle C1–C2–C3–C4 and specify 160° to set the first starting shape. Select the minimizer and allow it to run until it stops. Did it find the *anti* conformation (180°)? Record the energy. Repeat the process, starting with torsion angles of 0°, 45°, and 120° for the butane skeleton. Record the strain energies and report the final conformations that are formed in each case. What are your conclusions? Do your final results agree with those in the table?

If your minimizer rotated the two-eclipsed conformations (0° and 120°) to their closest staggered minima, you may have to restrict the minimizer to a single iteration in order to calculate their energies. This restriction calculates a **single-point energy**, and the energy of the structure is not minimized. If necessary, calculate the single-point energies of the eclipsed conformations and record your results.

¹ If your program does not have this feature, you can approximate the angles specified by constructing your starting molecules on the screen in a Z-shape for one and in a U-shape for the other.

The lesson here is that you may have to try several starting points to find the correct structure for the lowest energy conformation of a molecule! Do not blindly accept your first result, but look at it with the skeptical eye of a practiced chemist and test it further.

Optional Exercise. Record the single-point energies for every 30° rotation, starting at 0° and ending at 360°. When these energies are plotted against their angle, the plot should resemble the rotational energy curve shown for butane in most organic textbooks.

19B EXPERIMENT 19B

Cyclohexane Chair and Boat Conformations

In this exercise, we will investigate the chair and boat conformations of cyclohexane. Many programs will have these stored on disk as templates or fragments. If they are available as templates or fragments, you will need only to add hydrogens to the template. The chair is not difficult to build if you construct your cyclohexane on the screen in such a way that it resembles a chair (that is, just as you might draw it on paper). This crude construct will usually minimize to a chair. The boat is more difficult to construct. When you draw a crude boat on the screen, it will minimize to a *twist* boat, instead of the desired symmetrical boat.

Before you construct any cyclohexanes, construct a propane molecule. Minimize it and measure the CH and CC bond lengths and the CCC bond angle. Record these values; you will use them for reference.

Now construct a cyclohexane chair and minimize it. Measure the CH and CC bond lengths and the CCC angle in the ring. Compare these values to those of propane. What do you conclude? Rotate the molecule so that you view it end-on, looking down two of the bonds simultaneously (as in a Newman projection). Are all the hydrogens staggered? Rotate the chair and look at it from a different end-on angle. Are all the hydrogens still staggered? The van der Waals radius of a hydrogen atom is 1.20 Ångstroms. Hydrogen atoms that are closer than 2.40 Ångstroms apart will “touch” each other and create steric strain. Are any of the hydrogens in the cyclohexane chair close enough to cause steric strain? What are your conclusions?

Now construct a cyclohexane boat (from a template) and do not minimize it.¹ Measure the CH and CC bond lengths and the CCC bond angles at both the peaks and the lower corner of the ring. Compare these values to those of propane. Rotate the molecule so that you view it end-on, looking down the two parallel bonds on the sides of the boat. Are the hydrogens eclipsed or staggered? Now measure the distances between the various hydrogens on the ring, including the bowsprit-flagpole hydrogens and the axial and equatorial hydrogens on the side of the ring. Are any of the hydrogens generating steric strain?

Now minimize the boat to a twist boat and repeat all of the measurements. Write all of your conclusions about chairs, boats, and twist boats in your report.

¹ A single-point energy may be obtained, if you desire.

19C EXPERIMENT 19C

Substituted Cyclohexane Rings (Critical Thinking Exercises)

These exercises are designed to have you discover some not so obvious principles. Any conclusions and explanations that are requested should be recorded in your notebook.

Dimethylcyclohexanes. Using a cyclohexane template, construct *cis(a,a)*-1,3-dimethylcyclohexane, *cis(e,e)*-1,3-dimethylcyclohexane, and *trans(a,e)*-1,3-dimethylcyclohexane, and measure their energies. In the diaxial isomer, measure the distance between the two methyl groups. What do you conclude? Explain the result.

Similar comparisons can be made for the *cis*- and *trans*-1,2-dimethylcyclohexanes and the *cis*- and *trans*-1,4-dimethylcyclohexanes.

***cis*-1,4-Di-*tert*-butylcyclohexane.** Using hand drawings of chairs and boats, predict the expected conformation of this molecule. Then, construct *cis(a,e)*-1,4-di-*tert*-butylcyclohexane in a chair conformation, minimize it, and record its energy. Next, construct *cis(e,e)*-1,4-di-*tert*-butylcyclohexane in a boat conformation, placing the *tert*-butyl groups in equatorial positions at the peaks (puckered carbon atoms). Minimize this conformation to a twist boat and record its energy. Should we always expect chair conformations to have lower energy than boat conformations? Explain. What conformation do you predict for the *trans* stereoisomer?

***trans*-1,2-Dichloro and dibromocyclohexanes.** Build a model of *trans(a,a)*-1,2-dichlorocyclohexane, minimize it, and record its energy. Build a model of *trans(e,e)*-1,2-dichlorocyclohexane, minimize it, and record its energy. What is your conclusion? Now predict the result for the same two conformations of *trans*-1,2-dibromocyclohexane. When you have made a prediction, go ahead and model the two dibromo isomers and record the energies. What did you find? Explain the result. Do you think the result would be the same in a highly polar solvent?

Now construct the *cis*-1,2-dichloro and dibromocyclohexanes and compare their energies. Once again, explain what you find.

19D EXPERIMENT 19D

cis- and *trans*-2-Butene

Heats of hydrogenation for the three isomers of butene are given in the following table. Construct both *cis*- and *trans*-2-butene, minimize them, and report their energies. Which of these isomers has the lowest energy? Can you determine why?

Compound	ΔH (kcal/mol)	ΔH (kJ/mol)
<i>trans</i> -2-butene	-27.6	-115
<i>cis</i> -2-butene	-28.6	-120
1-butene	-30.3	-126

Now construct and minimize 1-butene. Record its energy. Obviously, 1-butene does not fit with the hydrogenation data. Molecular mechanics works quite well for *cis*- and *trans*-2-butene because they are very similar isomers. Both are 1,2-disubstituted alkenes. However, 1-butene is a monosubstituted alkene, and direct comparison to the 2-butenes cannot be made. The differences in the stability of mono- and disubstituted alkenes require that factors other than those used in molecular mechanics be used. These factors are caused by electronic and resonance differences. The molecular orbitals of the methyl groups interact with the pi bonds of the disubstituted alkenes (hyperconjugation) and help to stabilize them. Two such groups (as in 2-butene) are better than one (as in 1-butene). Therefore, although the bond lengths and angles come out pretty well for 1-butene, the energy derived for 1-butene does not directly compare to the energies of the 2-butenes. Molecular mechanics does not include terms that allow these factors to be included; it is necessary to use either semiempirical or *ab initio* quantum mechanical methods, which are based on molecular orbitals.

Computational Chemistry—*ab Initio* and Semiempirical Methods

In an earlier essay (“Molecular Modeling and Molecular Mechanics” that precedes Experiment 19), the application of **molecular mechanics** to solve chemical problems was discussed. Molecular mechanics is very good at giving estimates of the bond lengths and angles in a molecule. It can find the best geometry or conformation of a molecule. However, it requires the application of **quantum mechanics** to find good estimates of the thermodynamic, spectroscopic, and electronic properties of a molecule. In this essay, we will discuss the application of quantum mechanics to organic molecules.

Quantum mechanics computer programs can calculate heats of formation and the energies of transition states. The shapes of orbitals can be displayed in three dimensions. Important properties can be mapped onto the surface of a molecule. With these programs, the chemist can visualize concepts and properties in a way that the mind cannot readily imagine. Often this visualization is the key to understanding or to solving a problem.

INTRODUCTION TO TERMS AND METHODS

For you to solve the electronic structure and energy of a molecule, quantum mechanics requires that you formulate a wavefunction Ψ (psi) that describes the distribution of all the electrons within the system. The nuclei are assumed to have relatively small motions and to be essentially fixed in their equilibrium positions (Born–Oppenheimer approximation). The average energy of the system is calculated by using the Schrödinger equations as

$$E = \int \Psi^* H \Psi d\tau / \int \Psi^* \Psi d\tau$$

where H , the Hamiltonian operator, is a multiterm function that evaluates all the potential energy contributions (electron–electron repulsions and nuclear–electron attractions) and the kinetic energy terms for each electron in the system.

Because we can never know the true wavefunction Ψ for the molecule, we must guess at the nature of this function. According to the **Variation Principle**, a cornerstone idea in quantum mechanics, we can continue to guess at this function forever and never reach the true energy of the system, which will always be lower than our best guess. Because of the Variation Principle, we can formulate an approximate wavefunction and then consistently vary it until we minimize the energy of the system (as calculated using the Schrödinger equation). When we reach the variational minimum, the resulting wavefunction is often a good approximation of the system we are studying. Of course, you can’t just make any guess and get good results. It has taken theoretical chemists quite a few years to learn how to formulate both wavefunctions and Hamiltonian operators that yield results that agree quite closely with experimental data. Today, however, most methods for performing these calculations have been

well established, and computational chemists have devised easy-to-use computer programs, which can be used by any chemist to calculate molecular wavefunctions.

Molecular quantum-mechanical calculations can be divided into two classes: *ab initio* (Latin: “from the beginning” or “from first principles”) and *semiempirical*.

1. ***Ab initio* calculations** use the fully correct Hamiltonian operator for the system and attempt a complete solution without using any experimental parameters.
2. **Semiempirical calculations** generally use a simplified Hamiltonian operator and incorporate experimental data or a set of parameters that can be adjusted to fit experimental data.

Ab initio calculations require a great deal of computer time and memory, because every term in the calculations is evaluated explicitly. Semiempirical calculations have more modest computer requirements, allowing the calculations to be completed in a shorter time and making it possible to treat larger molecules. Chemists generally use semiempirical methods whenever possible, but it is useful to understand both methods when solving a problem.

SOLVING THE SCHRÖDINGER EQUATION

The Hamiltonian. The exact form of the Hamiltonian operator, which is a collection of potential energy (electrostatic attraction and repulsion terms) and kinetic energy terms, is now standardized and need not concern us here. However, all the programs require the **Cartesian coordinates** (locations in three-dimensional space) of all the atoms and a **connectivity matrix** that specifies which atoms are bonded and how (single, double, triple, H-bond, and so on). In modern programs, the user draws or constructs the molecule on the computer screen, and the program automatically constructs the atomic-coordinate and connectivity matrices.

The Wavefunction. It is not necessary for the user to construct or guess at a trial wavefunction—the program will do this. However, it is important to understand how the wavefunctions are formulated because the user frequently has a choice of methods. The complete molecular wavefunction is made up of a determinant of molecular orbitals:

$$\Psi = \begin{vmatrix} \phi_1(1) & \phi_2(1) & \phi_3(1) & \dots & \phi_n(1) \\ \phi_1(2) & \phi_2(2) & \phi_3(2) & \dots & \phi_n(2) \\ \phi_1(n) & \phi_2(n) & \phi_3(n) & \dots & \phi_n(n) \end{vmatrix}$$

The molecular orbitals $\phi_i(n)$ must be formulated from some type of mathematical function. They are usually made up of a **linear combination of atomic orbitals** χ_j (LCAO) from each of the atoms that make up the molecule.

$$\phi_i(n) = \sum_j c_{ji} \chi_j = c_1 \chi_1 + c_2 \chi_2 + c_3 \chi_3 \dots$$

This combination includes all the orbitals in the *core* and the *valence shell* of each atom in the molecule. The complete set of orbitals χ_j is called the **basis set** for the calculation. When an *ab initio* calculation is performed, most programs require the user to choose the basis set.


BASIS-SET ORBITALS

It should be apparent that the most obvious basis set to use for an *ab initio* calculation is the set of hydrogen-like atomic orbitals $1s$, $2s$, $2p$, and so on that we are all familiar with from atomic structure and bonding theory. Unfortunately, these “actual” orbitals present computational difficulties because they have radial nodes

when they are associated with the higher shells of an atom. As a result, a more convenient set of functions was devised by Slater. These **Slater-type orbitals (STOs)** differ from the hydrogen-like orbitals in that they have no radial nodes, but they have the same angular terms and overall shape. More importantly, they give good results (those that agree with experimental data) when used in semiempirical and *ab initio* calculations.

Slater-Type Orbitals. The radial term of an STO is an exponential function with the form $R_n = r^{(n-1)}e^{[-(Z-s)r/n]}$, where Z is the nuclear charge of the atom, and s is a “screening constant” that reduces the nuclear charge Z that is “seen” by an electron. Slater formulated a set of rules to determine the values of s that are required to produce orbitals that agree in shape with the customary hydrogen-like orbitals.

Radial Expansion and Contraction. A problem with simple STOs is that they do not have the ability to vary their radial size. Today it is common to use two or more simpler STOs so that expansion and contraction of the orbitals can occur during the calculation. For instance, if we take two functions such as $R(r) = r e^{(-\zeta r)}$ with different values of ζ , the larger value of ζ gives an orbital more contracted around the nucleus (an inner STO), and the smaller value of ζ gives an orbital extended further out from the nucleus (an outer STO). By using these two functions in different combinations, any size STO can be generated.

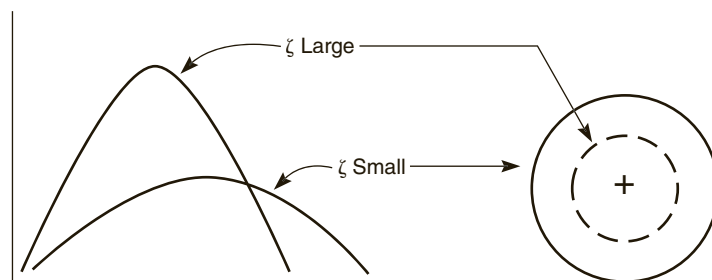


Variation of the radial size of an STO with the value of the exponent ζ (zeta).

Gaussian-Type Orbitals. The original Slater-type orbitals were eventually abandoned, and *simulated* STOs built from Gaussian functions were used. The most common basis set of this kind is the **STO-3G basis set**, which uses three Gaussian functions (3G) to simulate each one-electron orbital. A Gaussian function is of the type $R(r) = r e^{(-\alpha r^2)}$.

In the STO-3G basis set, the coefficients of the Gaussian functions are selected to give the best fit to the corresponding Slater-type orbitals. In this formulation, for instance, a hydrogen electron is represented by a single STO (a 1s type orbital) that is simulated by a combination of three Gaussian functions. An electron on any period 2 element (Li to Ne) will be represented by five STOs ($1s, 2s, 2p_x, 2p_y, 2p_z$), each simulated by three Gaussian functions. Each electron in a given molecule will have its own STO. (The molecule is literally built up by a series of one-electron orbitals. A spin function is also included so that no two of the one-electron orbitals are exactly the same.)

Split-Valence Basis Sets. A further step of evolution has made it now common to abandon attempts to simulate the hydrogen-like orbitals with STOs. Instead, an optimized combination of the Gaussian functions themselves is used for the basis set. The 3-21G basis set has largely replaced the STO-3G basis set for all but the largest molecules. The 3-21G symbolism means that three Gaussian functions are used for the wavefunction of each core electron, but the wavefunctions of the valence


electrons are “split” two-to-one (2:1) between inner and outer Gaussian functions, allowing the valence shell to expand or contract in size.

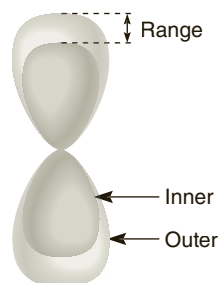

Split-valence orbitals.

A larger basis set (and one that requires more calculation time) is 6-31G, which uses six Gaussian “primitives” and a three-to-one split in the valence shell orbitals.

Polarization Basis Sets. Both the 3-21G and 6-31G basis sets can be extended to 3-21G* and 6-31G*. The star (*) indicates that these are **polarization sets**, in which the next higher type of orbital is included (for instance, a *p* orbital can be polarized by adding a *d* orbital function). Polarization allows deformation of the orbital toward the bond on one side of the atom.

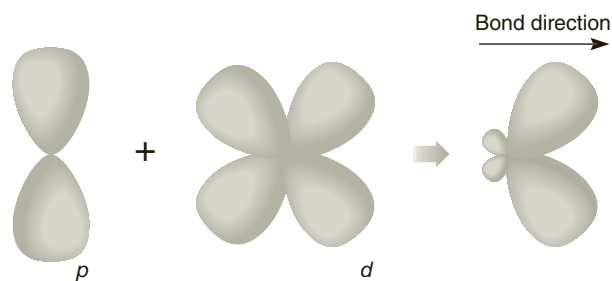

Polarization orbitals.

The largest basis set in current use is 6-311G*. Because it is computationally intensive, it is used only for **single-point calculations** (a calculation on a fixed geometry—no minimization performed). Other basis sets include the 6-31G** (which includes six *d* orbitals per atom instead of the usual five) and the 6-31+G* or 6-31++G* sets, which include diffuse *s* functions (electrons at a larger distance from the nucleus) to better deal with anions.

SEMIEMPIRICAL METHODS

It would be quite impossible to give a short and complete overview of the various semiempirical methods that have evolved over time. One must really get into the mathematical details of the method to understand what approximations have been made in each case and what kinds of empirical data have been included. In many of these methods, it is common to omit integrals that are expected (either from experience or for theoretical reasons) to have negligible values. Certain integrals are stored in a table and are not calculated each time the program is applied. For instance, the **frozen core approximation** is often used. This approximation assumes that the *completed shells* of the atom do not differ from one atom to another in the same period. All the core calculations are stored in a table, and they are simply looked up when needed. This makes the computation much easier to perform.

One of the more popular semiempirical methods in use today is AM-1. The parameters in this method work especially well for organic molecules. In fact, whenever possible, you should try to solve your problem using a semiempirical method such as AM-1 before you resort to an *ab initio* calculation. Also popular are MINDO/3 and MNDO, which are often found together in a computational package called MOPAC. If you are performing semiempirical calculations on inorganic molecules, you must make sure the method you use is optimized for transition metals. Two popular methods used by inorganic chemists wishing to involve metals in their calculations are PM-3 and ZINDO.

PICKING A BASIS SET FOR *AB INITIO* CALCULATIONS

When you perform an *ab initio* calculation, it is not always easy to know which basis set to use. Normally you should not use more complexity than is needed to answer your question or solve the problem. In fact, it may be desirable to determine the approximate geometry of the molecule using *molecular mechanics*. Many programs will allow you to use the result of a molecular mechanics **geometry optimization** as a starting point for an *ab initio* calculation. If possible, you should do so to save computational time.

Usually, 3-21G is a good starting point for an *ab initio* calculation, but if you have a very large molecule, you may wish to use STO-3G, a simpler basis set. Avoid doing geometry optimizations with the larger basis sets. Often you can do the geometry optimization first with 3-21G (or a semiempirical method) and then polish up the result with a **single-point energy** calculation with a larger basis set, such as 6-31G. You should “move up the ladder”: AM1 to STO-3G to 3-21G to 6-31G, and so on. If you don’t see any change in the results as you move up to successively more complex basis sets, it is generally fruitless to continue. If you include elements beyond period 2, use polarization sets (PM3 for semiempirical). Some programs have special sets for cations and anions or for radicals. If your result doesn’t match experimental results, you may not have used the correct basis set.

HEATS OF FORMATION

In classical thermodynamics, the **heat of formation**, ΔH_f , is defined as the energy consumed (endothermic reaction) or released (exothermic reaction) when a molecule is formed from its elements at standard conditions of pressure and temperature. The elements are assumed to be in their standard states.



Both *ab initio* and semiempirical programs calculate the energy of a molecule as its “heat of formation.” This heat of formation, however, is not identical to the thermodynamic function, and it is not always possible to make direct comparisons.

Heats of formation in semiempirical calculations are generally calculated in kcal/mole (1 kcal = 4.18 kJ) and are similar but not identical to the thermodynamic function. The AM1, PM3, and MNDO methods are parameterized by fitting them to a set of experimentally determined enthalpies. They are calculated from the binding energy of the system. The **binding energy** is the energy released when molecules are formed from their separated electrons and nuclei. The semiempirical heat of formation is calculated by subtracting atomic heats of formation from the binding energy. For most organic molecules, AM1 will calculate the heat of formation correctly to within a few kilocalories per mole.

In *ab initio* calculations, the heat of formation is given in **hartrees** (1 hartree = 627.5 kcal/mole = 2625 kJ/mole). In the *ab initio* calculation, the heat of formation is best defined as total energy. Like the binding energy, the **total energy** is the energy released when molecules are formed from their separated electrons and nuclei. This “heat of formation” always has a large negative value and does not relate well to the thermodynamic function.

Although these values do not relate directly to the thermodynamic values, they can be used to compare the energies of isomers (molecules of the same formula), such as *cis*- and *trans*-2-butene, or of tautomers, such as acetone in its enol and keto forms.


$$\Delta E = \Delta H_f(\text{isomer 2}) - \Delta H_f(\text{isomer 1})$$

It is also possible to compare the energies of balanced chemical equations by subtracting the energies of the products from the reactants.


$$\Delta E = [\Delta H_f(\text{product 1}) + \Delta H_f(\text{product 2})] - [\Delta H_f(\text{reactant 1}) + \Delta H_f(\text{reactant 2})]$$

GRAPHIC MODELS AND VISUALIZATION

Although the solution of the Schrödinger equation minimizes the *energy* of the system and gives a heat of formation, it also calculates the shapes and energies of all the molecular orbitals in the system. A big advantage of semiempirical and *ab initio* calculations, therefore, is the ability to determine the energies of the individual molecular orbitals and to plot their shapes in three dimensions. For chemists investigating chemical reactions, two molecular orbitals are of paramount interest: the HOMO and the LUMO.

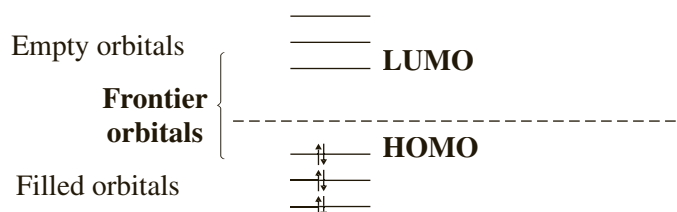

The **HOMO**, the highest occupied molecular orbital, is the last orbital in a molecule to be filled with electrons. The **LUMO**, the lowest unoccupied molecular orbital, is the first empty orbital in a molecule. These two orbitals are often called **frontier orbitals**.

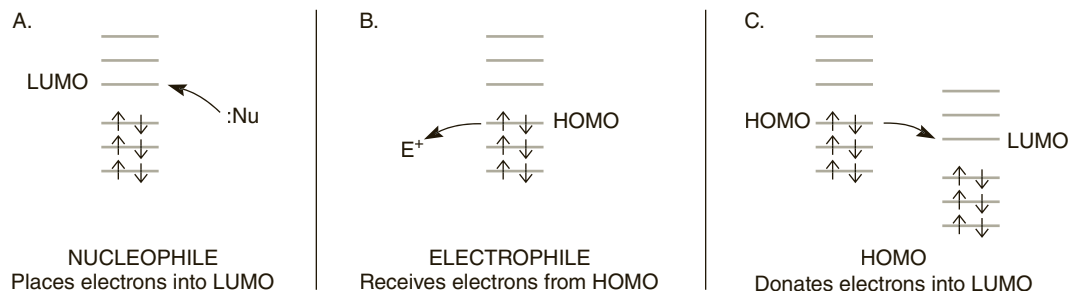


The frontier orbitals are similar to the valence shell of the molecule. They are where most of the chemical reactions occur. For instance, if a reagent is going to react with a Lewis base, the electron pair of the base must be placed into an empty

orbital of the acceptor molecule. The most available orbital is the LUMO. By examining the structure of the LUMO, one can determine the most likely spot where the addition will take place—usually at the atom where the LUMO has its biggest lobe. Conversely, if a Lewis acid attacks a molecule, it will bond to electrons that already exist in the molecule under attack. The most likely spot for this attack would be the atom where the HOMO has its biggest lobe (the electron density should be greatest at that site). Where it is not obvious which molecule is the electron pair donor, the HOMO that has the highest orbital energy will usually be the electron pair donor, placing electrons into the LUMO of the other molecule. The frontier orbitals, HOMO and LUMO, are where most chemical reactions occur.

SURFACES

Chemists use many kinds of handheld models to visualize molecules. A framework model best represents the angles, lengths, and directions of bonds. A molecule's size and shape are probably best represented by a space-filling model. In quantum mechanics, a model similar to the space-filling model can be generated by plotting a surface that represents all the points where the electron density of the molecule's wavefunction has a constant value. If this value is chosen correctly, the resulting surface will resemble the surface of a space-filling model. This type of surface is called an **electron-density surface**. The electron-density surface is useful for visualizing the size and shape of the molecule, but it does not reveal the position of the nuclei, bond lengths, or angles because you cannot see inside the surface. The electron-density value used to define this surface will be quite low because electron density falls off with increasing distance from the nucleus. If you choose a higher value of electron density when you plot this surface, a **bond-density surface** will be obtained. This surface will not give you an idea of the size or shape of the molecule, but it will reveal where the bonds are located, because the electron density will be higher where bonding is taking place.


Cyclopentane

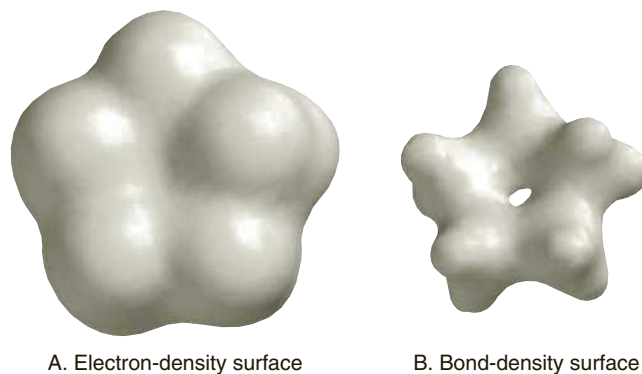

MAPPING PROPERTIES ONTO A DENSITY SURFACE

It is also possible to map a calculated property onto an electron-density surface. Because all three Cartesian coordinates are used to define the points on the surface, the property must be mapped in color, with the colors of the spectrum red–orange–yellow–green–blue representing a range of values. In effect, this is a four-dimensional plot ($x, y, z, +$ property mapped). One of the most common plots of this type is the **density–electrostatic potential**, or **density–elpot**, plot. The electrostatic potential is determined by placing a unit positive charge at each point

on the surface and measuring the interaction energy of this charge with the nuclei and electrons in the molecule. Depending on the magnitude of the interaction, that point on the surface is painted one of the colors of the spectrum. In the Spartan program, areas of high electron density are painted red or orange, and areas of lower electron density are plotted blue or green. When you view such a plot, the polarity of the molecule is immediately apparent.

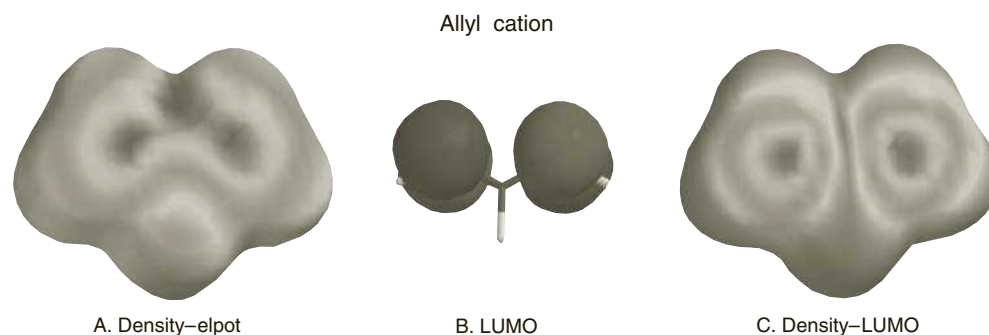

The second common type of mapping plots values of one of the frontier orbitals (either the HOMO or the LUMO) in color on the density surface. The color values plotted correspond to the value of the orbital where it intersects the surface. For a density-LUMO plot, for instance, the “hot spot” would be where the LUMO has its largest lobe. Because the LUMO is empty, this would be a bright blue area. In a density-HOMO plot, a bright red area would be the “hot spot.”

REFERENCES

Introductory

- Hehre, W. J.; Burke, L. D.; Shusterman, A. J.; Pietro, W. J. *Experiments in Computational Organic Chemistry*; Wavefunction, Inc.: Irvine, CA, 1993.
- Hehre, W. J.; Shusterman, A. J.; Nelson, J. E. *The Molecular Modeling Workbook for Organic Chemistry*; Wavefunction, Inc.: Irvine, CA 1998.
- Hypercube, Inc. *HyperChem Computational Chemistry*; HyperCube, Inc.: Waterloo, Ontario, Canada, 2014. The current version of this software is Version 8. Check the Hypercube website (www.hyper.com) for manuals and the most recent versions of the software. A student version of the software is also available.
- Shusterman, G. P.; Shusterman, A. J. Teaching Chemistry with Electron Density Models. *J. Chem. Educ.* **1997**, 74 (Jul), 771.
- Wavefunction, Inc. *PC-Spartan—Tutorial and User's Guide*; Wavefunction, Inc.: Irvine, CA, 2014. The current version of this software is Version 8. Check the Hypercube website (www.hyper.com) for manuals and the most recent versions of the software. A student version of the software is also available.

Advanced

- Clark, T. *Computational Chemistry*; Wiley-Interscience: New York, 1985.
- Fleming, I. *Frontier Orbitals and Organic Chemical Reactions*; John Wiley & Sons: New York, 1976.
- Fukui, K. Recognition of Stereochemical Paths by Orbital Interaction. *Accounts Chem. Res.* **1971**, 4, 57.
- Hehre, W. J.; Random, L.; Schleyer, P. v. R.; Pople, J. A. *Ab Initio Molecular Orbital Theory*; Wiley-Interscience: New York, 1986.
- Woodward, R. B.; Hoffmann, R. The Conservation of Orbital Symmetry. *Accounts Chem. Res.* **1968**, 1, 17.
- Woodward, R. B.; Hoffmann, R. *The Conservation of Orbital Symmetry*; Verlag Chemie: Weinheim, 1970.

Computational Chemistry

Semiempirical methods

Heats of formation

Mapped surfaces

REQUIRED READING

Review: The sections of your lecture textbook dealing with

20A: Alkene Isomers, Tautomerism, and Regioselectivity—the Zaitsev and Markovnikoff Rules

20B: Nucleophilic Substitution—Relative Rates of Substrates in S_N1 Reactions

20C: Acids and Bases—Inductive Effects

20D: Carbocation Stability

20E: Carbonyl Additions—Frontier Molecular Orbitals

New: Essay Computational Chemistry—*ab Initio* and Semiempirical Methods

SPECIAL INSTRUCTIONS

To perform this experiment, you must use computer software that can perform semiempirical molecular orbital calculations at the AM1 or MNDO level. In addition, the later experiments require a program that can display orbital shapes and map various properties onto an electron-density surface. Either your instructor will provide direction for using the software or you will be given a handout with instructions.

NOTES TO THE INSTRUCTOR

This series of computational experiments was devised using the programs PC Spartan and MacSpartan; however, it should be possible to use many other implementations of semiempirical molecular orbital theory. Some of the other capable programs for the PC and the Macintosh include HyperChem Release 5 and CAChe Workstation. You will need to provide your students with an introduction to your specific implementation. The introduction should show students how to build a molecule, how to select and submit calculations and surface models, and how to load and save files.

It is not intended that all these experiments be performed in a single session. They are intended to illustrate what you can do with computational chemistry but are not comprehensive. You may wish either to assign them with specific lecture

topics or to complement a particular experiment. Alternatively, you may wish to use them as patterns that students can use to devise their own computational procedures to solve a new problem.

For Experiments 20A and 20B, if your software will perform both AM1 (or a similar MNDO procedure) and calculations that include the effect of aqueous solvation (such as AM1-SM2), it may be instructive to have the students work in pairs. One student can perform gas-phase calculations, and the other can perform the same calculations, including the solvent effect. They can then compare results in their reports.

20A EXPERIMENT 20A


Heats of Formation: Isomerism, Tautomerism, and Regioselectivity

Part A. Isomerism

The stability of isomers may be directly compared by examining their heats of formation. In separate calculations, build models of *cis*-2-butene, *trans*-2-butene, and 1-butene. Submit each of these to AM1 calculation of the energy (heat of formation). Use the geometry optimization option in each case to find the best possible energy for each isomer. What do your results suggest? Do they agree with the experimental data given in Experiment 19D?

Part B. Acetone and its Enol

In this exercise, we will compare the energies of a pair of tautomers using the heats of formation calculated by the semiempirical AM1 method. These two tautomers can be directly compared because they have the same molecular formula: C₃H₆O. Most organic textbooks discuss the relative stability of ketones and their tautomeric enol forms. For acetone, there are two tautomers in equilibrium:

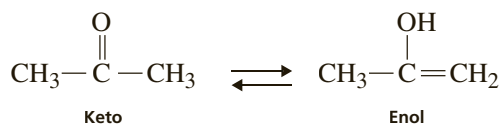


In separate calculations, build models of both acetone and its enol. Submit each model to AM1 calculation of the energy (heat of formation). Use the geometry optimization option in each case to find the best possible energy for each tautomer.

Experimental results indicate that there is very little enol (<0.0002%) in equilibrium with acetone. Do your calculations suggest a reason?


Part C. Regioselectivity

Ionic addition reactions of alkenes are quite regioselective. For instance, adding concentrated HCl to 2-methylpropene produces largely 2-chloro-2-methylpropane and a much smaller amount of 1-chloro-2-methylpropane. This can be explained by examining the energies of the two carbocation A and B that can be formed by adding a proton in the first step of the reaction:

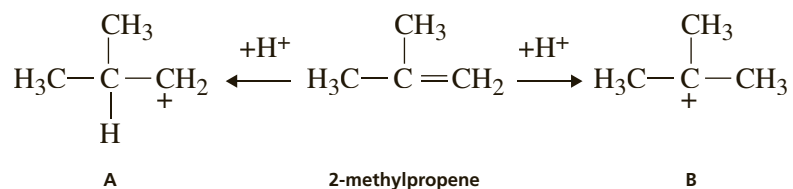

This first step (adding a proton) is the rate-determining step of the reaction, and it is expected that the activation energies for forming these two intermediates will reflect their relative energies. That is, the activation energy leading to the lower-energy intermediate will be lower than the activation energy leading to the intermediate that has higher energy. Because of this energy difference, the reaction will predominantly follow the pathway that passes through the lower-energy intermediate. Because the two carbocations are isomers and because both are formed from the same starting material, a direct comparison of their energies (heats of formation) will determine the main course of the reaction.

In separate calculations, build models of the two carbocations and submit them to AM1 calculations of their energies. Use a geometry optimization. When you build the models, most programs will require you to build the skeleton of the hydrocarbon that is closest in structure to the carbocation and then to delete the required hydrogen *and its free valence*.

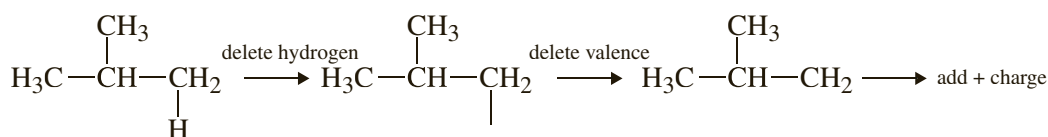


Remember also to assign a positive charge to the molecule before submitting it to calculation. This is usually done in the menus where you select the type of calculation. Compare your results for the two calculations. Which carbocation will lead to the major product? Do your results agree with the prediction made by Markovnikoff's Rule?

20B EXPERIMENT 20B

Heats of Reaction: S_N1 Reaction Rates

In this experiment, we will attempt to determine the relative rates of selected substrates in the S_N1 reaction. The effect of the degree of substitution will be examined for the following compounds:

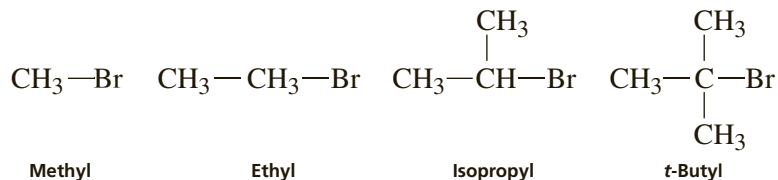

Because the four carbocations are not isomers, we cannot compare their heats of formation directly. To determine the relative rates at which these compounds react, we must determine the *activation energy* required to form the carbocation intermediate in each case. Ionization is the rate-determining step, and we will assume that the activation energy for each ionization should be *similar in magnitude* (Hammond Postulate) to the calculated energy difference between the alkyl halide and the two ions that it forms.



$$\Delta E_{\text{activation}} \cong \Delta H_{\text{f}}(\text{products}) - \Delta H_{\text{f}}(\text{reactants}) \quad [2]$$

$$\Delta E_{\text{activation}} \cong \Delta H_{\text{f}}(\text{R}^+) + \Delta H_{\text{f}}(\text{Br}^-) - \Delta H_{\text{f}}(\text{RBr}) \quad [3]$$

Because the energy of the bromide ion is a constant, it could be omitted from the calculation, but we will include it because it must be computed only once.

Part A. Ionization Energies

Using the AM1 semiempirical level of calculation, compute the energies (heats of formation) of each of the starting materials and record them. Next, compute the energies of each of the carbocations that would result from the ionization of each substrate—follow the instructions given in Part C of Experiment 20A—and record the results. Be sure to add the positive charge. Finally, compute the energy of the bromide ion, remembering to delete the free valence and add a negative charge. Once all the calculations have been performed, use Equation 3 to calculate the energy required to form the carbocation in each case. What do you conclude about the relative rates of the four compounds?


Part B. Solvation Effects (Optional)

The calculations you performed in Part A did not take the effect of solvation of the ions into account. At your instructor's option (and if you have the correct software), you may be required to repeat your calculations using a computational method that includes stabilization of the ions by solvation. Will solvation increase or decrease the ionization energies? Which will be solvated more, the reactants or the products of the ionization step? What do you conclude from your results?

20C EXPERIMENT 20C

Density–Electrostatic Potential Maps: Acidities of Carboxylic Acids

In this experiment, we will compare the acidities of acetic, chloroacetic, and trichloroacetic acid. This experiment could be approached in the same fashion as the relative rates in Experiment 20B, using the ionization energies to determine the relative acidities.

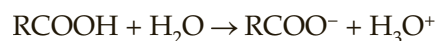

$$\Delta E = [\Delta H_{\text{f}}(\text{RCOO}^-) + \Delta H_{\text{f}}(\text{H}_3\text{O}^+)] - [\Delta H_{\text{f}}(\text{RCOOH}) + \Delta H_{\text{f}}(\text{H}_2\text{O})]$$

In fact, the water and hydronium ion terms could be omitted because they would be constant in each case.

Instead of calculating the ionization energies, we will use a more visual approach involving a property map. Set up an AM1 geometry optimization calculation for each of the acids. In addition, request that an electron-density surface be calculated with the electrostatic potential mapped onto this surface in color. In this procedure, the program plots the density surface and determines the electron density at each point by placing a test positive charge there and determining the coulomb interaction. The surface is colored using the colors of the spectrum—blue is used for positive areas (low electron density), and red is used for more negative areas (high electron density). This plot will show the polarization of the molecule.

When you have finished the calculations, display all three maps on the screen at the same time. To compare them, you must adjust them all to the same set of color values. This can be done by observing the maximum and minimum values for each map in the surface display menus. Once you have all six values (save them), determine which two numbers give you the maximum and minimum values. Return to the surface plot menu for each of the molecules and readjust the limits of the color values to the same maximum and minimum values. Now the plots will all be adjusted to identical color scales. What do you observe for the carboxyl protons of acetic acid, chloroacetic acid, and trichloroacetic acid? The three minimum values that you saved can be compared to determine the relative electron density at each proton.

20D EXPERIMENT 20D

Density–Electrostatic Potential Maps: Carbocations

Part A. Increasing Substitution

In this experiment, we will use a density map to determine how well a series of carbocations disperses the positive charge. According to theory, increasing the number of alkyl groups attached to the carbocation center helps to spread out the charge (through hyperconjugation) and lowers the energy of the carbocation. We approached this problem from a computational (numerical) angle in Experiment 20B. Now we will prepare a visual solution to the problem.

Begin by performing an AM1 geometry optimization on methyl, ethyl, isopropyl, and *tert*-butyl carbocations. These carbocations are built as described in Part C of Experiment 20A. Don't forget to specify that each one has a positive charge. Also select a density surface for each one with the electrostatic potential mapped onto the surface.

When the calculations are completed, display all four density–electrostatic potential maps on the same screen and adjust the color values to the same range as described in Experiment 20C. What do you observe? Is the positive charge as localized in the *tert*-butyl carbocation as in its methyl counterpart?

Part B. Resonance

Repeat the computational experiment described in Part A, using density–electrostatic potential maps for the allyl and benzyl carbocations. These two experiments can be performed without displaying them both on the same screen. What do you observe about the charge distribution in these two carbocations?


20E EXPERIMENT 20E

Density–LUMO Maps: Reactivities of Carbonyl Groups

In this experiment, we will investigate how frontier molecular orbital theory applies to the reactivity of a carbonyl compound. Consider the reaction of a nucleophile such as hydride or cyanide with a carbonyl compound.

According to frontier molecular orbital theory (see the section “Graphic Models and Visualization” in the essay that precedes this experiment), the nucleophile,


which is donating electrons, must place them in an empty orbital of the carbonyl. Logically, this empty orbital would be the LUMO—the Lowest (energy) Unoccupied Molecular Orbital.

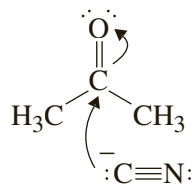

Make a model of acetone and submit it to an AM1 calculation with geometry optimization. Also select two surfaces to display, the LUMO and a mapping of the LUMO on a density surface.

When the calculations are finished, display both surfaces on the screen at the same time. Where is the biggest lobe of the LUMO, on carbon or on oxygen? Where does the nucleophile attack? The density–LUMO surface displays the same thing, but with color coding. This plot shows a blue spot on the surface where the LUMO has its greatest density (largest lobe).

Next, continue this experiment by calculating the LUMO and the density–LUMO plots for the ketones 2-cyclohexenone and norbornanone.

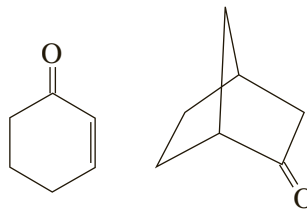

Where are the reactive sites in cyclohexenone? According to the literature, strong bases, such as Grignard reagents, attack the carbonyl, and weaker bases or better nucleophiles, such as amines, attack the beta carbon of the double bond, performing a conjugate addition. Can you explain this? Will a nucleophile attack norbornanone from the exo (top) or the endo (bottom) face of the molecule? See Experiment 33 for the answer.

Properties and Reactions of Organic Compounds

Reactivities of Alkyl Halides

S_N1/S_N2 reactions


Relative rates

Reactivities

The reactivities of alkyl halides in nucleophilic substitution reactions depend on two important factors: reaction conditions and substrate structure. The reactivities of several substrate types will be examined under both S_N1 and S_N2 reaction conditions in this experiment.



Sodium Iodide or Potassium Iodide in Acetone

A reagent composed of sodium iodide or potassium iodide dissolved in acetone is useful in classifying alkyl halides according to their reactivity in an S_N2 reaction. Iodide ion is an excellent nucleophile, and acetone is a nonpolar solvent. The tendency to form a precipitate increases the completeness of the reaction. Sodium iodide and potassium iodide are soluble in acetone, but the corresponding bromides and chlorides are not soluble. Consequently, as bromide ion or chloride ion is produced, the ion is precipitated from the solution. According to LeChâtelier's Principle, the precipitation of a product from the reaction solution drives the equilibrium toward the right; such is the case in the reaction described here:

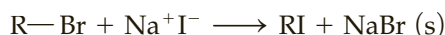

Silver Nitrate in Ethanol

A reagent composed of silver nitrate dissolved in ethanol is useful in classifying alkyl halides according to their reactivity in an S_N1 reaction. Nitrate ion is a poor nucleophile, and ethanol is a moderately powerful ionizing solvent. The silver ion, because of its ability to coordinate the leaving halide ion to form a silver halide precipitate, greatly assists the ionization of the alkyl halide. Again, a precipitate as one of the reaction products also enhances the reaction.

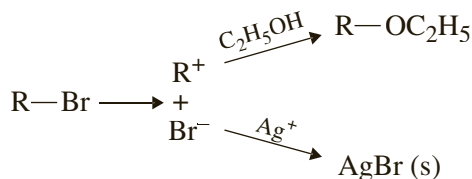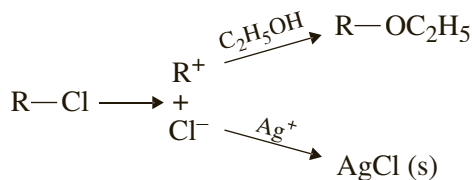

REQUIRED READING

Before beginning this experiment, review the chapters dealing with nucleophilic substitution in your lecture textbook.

SPECIAL INSTRUCTIONS

Some compounds used in this experiment, particularly crotyl chloride and benzyl chloride, are powerful lachrymators. **Lachrymators** cause eye irritation and the formation of tears.

CAUTION



Because some of these compounds are lachrymators, perform these tests in a hood. Be careful to dispose of the test solutions in a waste container marked for halogenated organic waste. After testing, rinse the test tubes with acetone and pour the contents into the same waste container.

SUGGESTED WASTE DISPOSAL

Dispose of all the halide wastes into the container marked for halogenated waste. Any acetone washings should also be placed in the same container.

NOTES TO THE INSTRUCTOR

Each of the halides should be checked with NaI/acetone and AgNO₃/ethanol to test for their purity before the class performs this experiment. If molecular modeling software is available, you may wish to assign the exercises included at the end of this experiment.

An alternative approach¹ for conducting this experiment is to restrict the list of test compounds to the following five substrates: 1-chlorobutane, 1-bromobutane, 2-chlorobutane, 2-bromobutane, and 2-chloro-2-methylpropane (*tert*-butyl chloride). If conducted in this way, one can simplify the experiment by eliminating the allylic, benzylic, and halocycloalkanes. This experiment can best be used if assigned *before* the S_N1 and S_N2 reactions have been discussed in lecture! An excellent and meaningful guided-inquiry experience can then be achieved by having students submit their results to a campus discussion board, such as BlackBoard or Canvas, prior to any discussion of the results by the instructor. Once the class results have been posted on BLACKBOARD or CANVAS, have the students study the class data to look for patterns. Encourage the class to try to “discover” how the reactivities in the sodium iodide/acetone and silver nitrate/ethanol depend on the substrate structure and the leaving group.

¹ This approach was suggested and utilized successfully by Professor Emily Borda, Department of Chemistry, Western Washington University, Bellingham, WA 98225. The authors wish to thank Professor Borda for her excellent contribution.

PROCEDURE

Part A. Sodium Iodide in Acetone

The Experiment

Label a series of ten clean and dry test tubes (10 × 75 mm test tubes may be used) from 1 to 10. In each test tube, place 2 mL of a 15% NaI-in-acetone solution. Now add 4 drops of one of the following halides to the appropriate test tube: (1) 2-chlorobutane, (2) 2-bromobutane, (3) 1-chlorobutane, (4) 1-bromobutane, (5) 2-chloro-2-methylpropane (*t*-butyl chloride), (6) crotyl chloride $\text{CH}_3\text{CH}=\text{CHCH}_2\text{Cl}$ (see Special Instructions), (7) benzyl chloride (α -chlorotoluene) (see Special Instructions), (8) bromobenzene, (9) bromocyclohexane, and (10) bromocyclopentane. Make certain you return the dropper to the proper container to avoid cross-contaminating these halides.

Reaction at Room Temperature

After adding the halide, shake the test tube² well to ensure adequate mixing of the alkyl halide and the solvent. Record the times needed for any precipitate or cloudiness to form.

Reaction at Elevated Temperature

After about 5 minutes, place any test tubes that do not contain a precipitate in a 50°C water bath. Be careful not to allow the temperature of the water bath to exceed 50°C, because the acetone will evaporate or boil out of the test tube. After about 1 minute of heating, cool the test tubes to room temperature and note whether a reaction has occurred. Record the results.

Observations

Generally, reactive halides give a precipitate within 3 minutes at room temperature, moderately reactive halides give a precipitate when heated, and unreactive halides do not give a precipitate, even after being heated. Ignore any color changes.

Report

Record your results in tabular form in your notebook. Explain why each compound has the reactivity you observed. Explain the reactivities in terms of structure. Compare relative reactivities for compounds of similar structure.

Part B. Silver Nitrate in Ethanol

The Experiment

Label a series of ten clean and dry test tubes from 1 to 10, as described in the previous section. Add 2 mL of a 1% ethanolic silver nitrate solution to each test tube. Now add 4 drops of the appropriate halide to each test tube, using the same numbering scheme indicated for the sodium iodide test. To avoid cross-contaminating these halides, return the dropper to the proper container.

Reaction at Room Temperature

After adding the halide, shake the test tube well to ensure adequate mixing of the alkyl halide and the solvent. After thoroughly mixing the samples, record the times needed for any precipitate or cloudiness to form. Record your results as dense precipitate, cloudiness, or no precipitate/cloudiness.

² Do not use your thumb or a stopper. Instead, hold the top of the test tube between the thumb and index finger of one hand and “flick” the bottom of the test tube using the index finger of your other hand.

Reaction at Elevated Temperature

After about 5 minutes, place any test tubes that do not contain a precipitate or cloudiness in a hot water bath at about 100°C. After about 1 minute of heating, cool the test tubes to room temperature and note whether a reaction has occurred. Record your results as dense precipitate, cloudiness, or no precipitate/cloudiness.

Observations

Reactive halides give a precipitate (or cloudiness) within 3 minutes at room temperature, moderately reactive halides give a precipitate (or cloudiness) when heated, and unreactive halides do not give a precipitate, even after being heated. Ignore any color changes.

Report

Record your results in tabular form in your notebook. Explain why each compound has the reactivity that you observed. Explain the reactivities in terms of structure. Compare relative reactivities for compounds of similar structure.

MOLECULAR MODELING (OPTIONAL)

Many points developed in this experiment can be confirmed through the use of molecular modeling. The following experiments were developed with PC Spartan. It should be possible to use other software, but the instructor may have to make some modifications.

S_N1 Reactivities

Part One. The rate of an S_N1 reaction is related to the energy of the carbocation intermediate that is formed in the rate-determining ionization step of the reaction. It is expected that the activation energy required to form an intermediate is close to the energy of the intermediate. When two intermediates are compared, the activation energy leading to the intermediate of lower energy is expected to be lower than the activation energy leading to the intermediate of higher energy. The easier it is to form the carbocation, the faster the reaction will proceed. An AM1 semiempirical method for determining the approximate energies of carbocation intermediates is described in Experiment 20B. Complete the computational exercises in Experiment 20B, and compare the calculated results to the experimental results you obtained in this experiment. Do the experimental results parallel the calculated results?

Part Two. Using the density–elpot surface plot described in Experiment 20D, it is possible to compare the amount of charge delocalization in various carbocations through a visualization of the ions. Complete Experiment 20D, and determine whether the charge distributions (delocalization) are what you would expect for the series of carbocations studied.

Part Three. The benzyl (and allyl) halides are a special case; they have resonance. To see how the charge is delocalized in the benzyl carbocation, request two plots: the electrostatic potential mapped onto a density surface and the LUMO mapped onto a density surface. Submit these for calculation at the AM1 semiempirical level. On a piece of paper, draw the resonance-contributing structures for the benzyl cation. Do the computational results agree with the conclusions you draw from your resonance hybrid?


Part Four. Repeat the calculation outlined in Part Three for the benzyl cation; however, in this calculation, turn the CH₂ group so that its hydrogens are perpendicular to the plane of the benzene ring. Compare your results to those obtained in Part Three.

S_N2 Reactivities

The problem in the S_N2 reaction is not an electric one, but rather a steric problem. Using the AM1 semiempirical method, request a LUMO surface and a density surface for each substrate. The simplest way to visualize the steric problem is to plot the LUMO inside a density surface mapped as a net or a transparent surface. Now imagine having to attack the back lobe of the LUMO. Compare bromomethane, 2-bromo-2-methylpropane (*tert*-butyl bromide), and 1-bromo-2,2-dimethylpropane (neopentyl bromide). Is there any electron density (atoms) in the way of the nucleophile? Request and calculate another surface, mapping the LUMO onto the density surface. What are your conclusions? Can you find the “hot spot” where the nucleophile will attack? Is there any steric hindrance?

QUESTIONS

1. In the tests with sodium iodide in acetone and silver nitrate in ethanol, why should 2-bromobutane react faster than 2-chlorobutane?
2. Why is benzyl chloride reactive in both tests, whereas bromobenzene is unreactive?
3. When benzyl chloride is treated with sodium iodide in acetone, it reacts much faster than 1-chlorobutane, even though both compounds are primary alkyl chlorides. Explain this rate difference.
4. 2-Chlorobutane reacts much more slowly than 2-chloro-2-methylpropane in the silver nitrate test. Explain this difference in reactivity.
5. Bromocyclopentane is more reactive than bromocyclohexane when heated with sodium iodide in acetone. Explain this difference in reactivity.
6. How do you expect the following series of compounds to compare in behavior in the two tests?

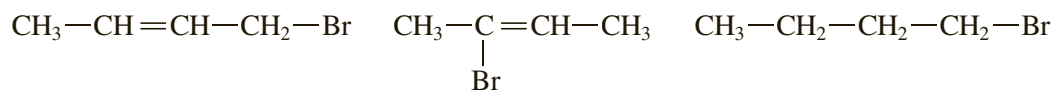

Nucleophilic Substitution Reactions: Competing Nucleophiles

Nucleophilic substitution

Heating under reflux

Extraction

Gas chromatography



NMR spectroscopy

In this experiment, you will compare the relative nucleophilicities of chloride ions and bromide ions toward each of the following alcohols: 1-butanol (*n*-butyl alcohol), 2-butanol (*sec*-butyl alcohol), and 2-methyl-2-propanol (*t*-butyl alcohol). The two nucleophiles will be present at the same time in each reaction, in equimolar concentrations, and they will be competing for substrate. A protic solvent is used in these reactions.

In general, alcohols do not react readily in simple nucleophilic displacement reactions. If they are attacked by nucleophiles directly, hydroxide ion, a strong base, must be displaced. Such a displacement is not energetically favorable and cannot occur to any reasonable extent:



To avoid this problem, you must carry out nucleophilic displacement reactions on alcohols in acidic media. In a rapid initial step, the alcohol is protonated; then water, a stable molecule, is displaced. This displacement is energetically favorable, and the reaction proceeds in high yield:

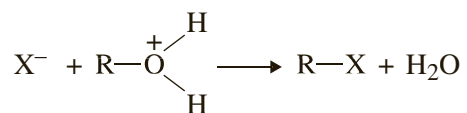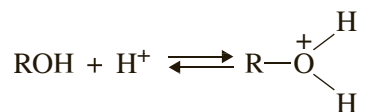

Once the alcohol is protonated, it reacts by either the S_N1 or the S_N2 mechanism, depending on the structure of the alkyl group of the alcohol. For a brief review of these mechanisms, consult the chapters on nucleophilic substitution in your lecture textbook.

You will analyze the products of the three reactions in this experiment by a variety of techniques to determine the relative amounts of alkyl chloride and alkyl bromide formed in each reaction. That is, using equimolar concentrations of

chloride ions and bromide ions reacting with 1-butanol, 2-butanol, and 2-methyl-2-propanol, you will determine which ion is the better nucleophile. In addition, you will determine for which of the three substrates (reactions) this difference is important and whether an S_N1 or S_N2 mechanism predominates in each case.

REQUIRED READING

Review: Techniques 1 through 6

Technique 7 Reaction Methods, Section 7.2, 7.4, and 7.8

Technique 12 Extractions, Separations, and Drying Agents, Sections 12.5, 12.9, and 12.11

Technique 21 Nuclear Magnetic Resonance Spectroscopy

Technique 22 Gas Chromatography

Before beginning this experiment, review the appropriate chapters on nucleophilic substitution in your lecture textbook.

SPECIAL INSTRUCTIONS

Each student will carry out the reaction with 2-methyl-2-propanol. Your instructor will also assign you either 1-butanol or 2-butanol. By sharing your results with other students, you will be able to collect data for all three alcohols. You should begin this experiment with Experiment 22A. During the lengthy reflux period, you will be instructed to go on to Experiment 22B. When you have prepared the product of that experiment, you will return to complete Experiment 22A. To analyze the results of both experiments, your instructor will assign specific analysis procedures in Experiment 22C that the class will accomplish.

The solvent–nucleophile medium contains a high concentration of sulfuric acid. Sulfuric acid is corrosive; be careful when handling it.

In each experiment, the longer your product remains in contact with water or aqueous sodium bicarbonate, the greater the risk that your product will decompose, leading to errors in your analytical results. Before coming to class, prepare so that you know exactly what you are supposed to do during the purification stage of the experiment.

SUGGESTED WASTE DISPOSAL

When you have completed the three experiments and all the analyses have been completed, discard any remaining alkyl halide mixture in the organic waste container marked for the disposal of halogenated substances. All aqueous solutions produced in this experiment should be disposed of in the container for aqueous waste.

NOTES TO THE INSTRUCTOR

The solvent–nucleophile medium must be prepared in advance for the entire class. Use the following procedure to prepare the medium.

This procedure will provide enough solvent–nucleophile medium for about 10 students (assuming no spillage or other types of waste). Place 100 g of ice in a 500-mL Erlenmeyer flask and carefully add 76 mL concentrated sulfuric acid. Carefully weigh 19.0 g ammonium chloride and 35.0 g ammonium bromide into a beaker. Crush any lumps of the reagents to powder and then, using a powder funnel, transfer the halides to an Erlenmeyer flask. Carefully add the sulfuric acid mixture to the ammonium salts a little at a time. Swirl the mixture vigorously to dissolve the salts. It will probably be necessary to heat the mixture on a steam bath or a hot plate to achieve total solution. Keep a thermometer in the mixture and make sure that the temperature does not exceed 45°C. If necessary, you may add as much as 10 mL of water at this stage. Do not worry if a few small granules do not dissolve. When solution has been achieved, pour the solution into a container that can be kept warm until all students have taken their portions. The temperature of the mixture must be maintained at about 45°C to prevent precipitation of the salts. Be careful that the solution temperature does not exceed 45°C, however. Place a 10-mL or 20-mL calibrated pipette fitted with a pipette pump in the mixture. The pipette is always left in the mixture to keep it warm.

Be certain that the *tert*-butyl alcohol has been melted before the beginning of the laboratory period.

The gas chromatograph should be prepared as follows: column temperature, 100°C; injection and detector temperature, 130°C; carrier gas flow rate, 50 mL/min. The recommended column is 8 feet long, with a stationary phase such as Carbowax 20M. If you wish to analyze the products from the reaction of *tert*-butyl alcohol (Exp. 22B) by gas chromatography, be sure that the *tert*-butyl halides do not undergo decomposition under the conditions set for the gas chromatograph. *tert*-Butyl bromide is susceptible to elimination. It will be necessary to determine retention times of the pure compounds to determine the order of elution of the halides.

Unless the samples are analyzed by gas chromatography immediately after preparing them, it is essential that the samples be stored in leak-proof vials. We have found GC-MS vials to be ideal for this purpose.


22A EXPERIMENT 22A

Competing Nucleophiles with 1-Butanol or 2-Butanol

PROCEDURE

Apparatus

Assemble an apparatus for reflux using a 20-mL or 25-mL round-bottom flask, a reflux condenser, and a drying tube, as shown in the figure. Loosely insert dry glass wool into the drying tube and then add water dropwise onto the glass wool until it is partially moistened. The moistened glass wool will trap the hydrogen chloride and hydrogen bromide gases produced during the reaction. As an alternative, you can use an external gas trap as described in Technique 7, Section 7.7, Part B. Do not place the round-bottom flask into the aluminum block until the reaction mixture has been added to the flask. Six Pasteur pipettes, two 3-mL conical vials with Teflon cap liners, and a 5-mL conical vial with a Teflon liner should also be assembled for use. All pipettes and vials should be clean and dry.

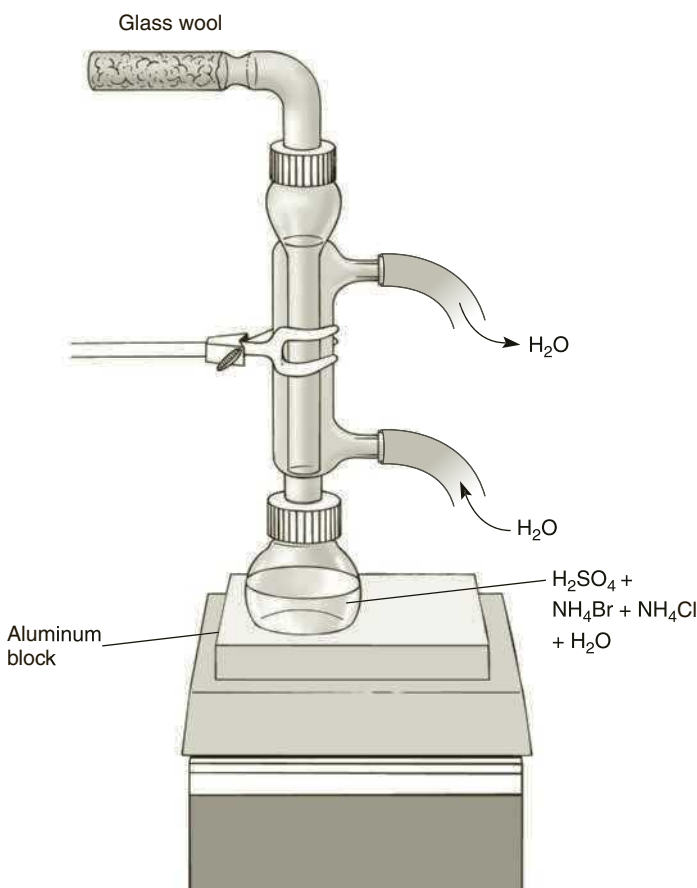

Apparatus for reflux.

CAUTION

The solvent–nucleophile medium contains a high concentration of sulfuric acid. This liquid will cause severe burns if it touches your skin.

Preparation of Reagents

If a calibrated pipette fitted with a pipette pump is provided, you may adjust the pipette to 10 mL and deliver the solvent–nucleophile medium directly into your 20-mL round-bottom flask (temporarily placed in a beaker for stability). Alternatively, you may use a warm 10-mL graduated cylinder to obtain 10.0 mL of the solvent–nucleophile medium. The graduated cylinder must be warm to prevent precipitation of the salts. Heat it by running hot water over the outside of the cylinder or by putting it in the oven for a few minutes. Immediately pour the mixture into the round-bottom flask. With either method, a small portion of the salts in the flask may precipitate as the solution cools. Do not worry about this; the salts will redissolve during the reaction.

Reflux

Assemble the apparatus shown in the figure. Using the following procedure, add 0.75 mL of 1-butanol (*n*-butyl alcohol) or 0.75 mL of 2-butanol (*sec*-butyl alcohol), depending on which alcohol you were assigned, to the solvent–nucleophile mixture contained in the reflux apparatus. Dispense the alcohol from the automatic pipette or dispensing pump into a test tube. Remove the drying tube and, with a 9-inch Pasteur pipette, dispense the alcohol directly into the round-bottom flask by inserting the Pasteur pipette into the opening of the condenser. Also add an inert boiling stone.¹ Replace the drying tube and start circulating the cooling water. Lower the reflux apparatus so that the round-bottom flask is in the aluminum block, as shown in the figure. Adjust the heat so that this mixture maintains a *gentle* boiling action. For 1-butanol, the aluminum block temperature should be about 140°C, and for 2-butanol, the temperature should be about 120°C. Be careful to adjust the reflux ring, if one is visible, so that it remains in the lower fourth of the condenser. Violent boiling will cause loss of product. Continue heating the reaction mixture containing 1-butanol for 75 minutes. Heat the mixture containing 2-butanol for 60 minutes. During this heating period, go on to Experiment 22B and complete as much of it as possible before returning to this procedure.

Purification

When the period of reflux has been completed, discontinue heating, lift the apparatus out of the aluminum block, and allow the reaction mixture to cool. Do not remove the condenser until the flask is cool. Be careful not to shake the hot solution

¹ Do not use calcium carbonate-based stones or Boileezers, because they will partially dissolve in the highly acidic reaction mixture.

as you lift it from the heating block or a violent boiling and bubbling action will result; this could allow material to be lost out the top of the condenser. After the mixture has cooled for about 5 minutes, immerse the round-bottom flask (with condenser attached) in a beaker of cold tap water (no ice) and wait for this mixture to cool down to room temperature.

There should be an organic layer present at the top of the reaction mixture. Add 0.75 mL of pentane to the mixture and *gently* swirl the flask. The purpose of the pentane is to increase the volume of the organic layer so that the following operations are easier to accomplish. Using a Pasteur pipette, transfer most (about 7 mL) of the bottom (aqueous) layer to another container. Be careful that all of the top organic layer remains in the boiling flask. Transfer the remaining aqueous layer and the organic layer to a 3-mL conical vial, taking care to leave behind any solids that may have precipitated. Allow the phases to separate and remove the bottom (aqueous) layer using a Pasteur pipette.

NOTE: For the following sequence of steps, be certain to be well prepared. If you find that you are taking longer than 5 minutes to complete the entire extraction sequence, you probably have affected your results adversely!

Add 1.0 mL of water to the vial and *gently* shake this mixture. Allow the layers to separate and remove the aqueous layer, which is still on the bottom. Extract the organic layer with 1–2 mL of saturated sodium bicarbonate solution and remove the bottom aqueous layer.

Drying

Using a clean dry Pasteur pipette, transfer the remaining organic layer into a small test tube (10 × 75 mm) and dry over anhydrous granular sodium sulfate (see Technique 12, Section 12.9). Transfer the dry halide solution with a clean, dry Pasteur pipette to a small, dry leak-proof vial, taking care not to transfer any solid.² *Be sure the cap is screwed on tightly.* Do not store the liquid in a container with a cork or a rubber stopper, because these will absorb the halides. This sample can now be analyzed by as many of the methods in Experiment 22C as your instructor indicates. If possible, analyze the sample on the same day.

22B EXPERIMENT 22B

Competing Nucleophiles with 2-Methyl-2-Propanol

PROCEDURE

Place 6.0 mL of the solvent–nucleophile medium into a 15-mL centrifuge tube, using the same procedure described in the “Preparation of Reagents” section at the beginning of Experiment 22A. Place the centrifuge tube in cold tap water and wait until a few crystals of ammonium halide salts just begin to appear. Using an automatic pipette or dispensing pump, transfer 1.0 mL of 2-methyl-2-propanol (*tert*-butyl alcohol, mp 25°C) to the 15-mL centrifuge tube. Replace the cap and make sure that it doesn’t leak.

² We have found GC-MS vials ideal for this purpose.

CAUTION

The solvent–nucleophile mixture contains concentrated sulfuric acid.

Shake the tube vigorously, venting occasionally, for 5 minutes (use gloves). Any solids that were originally present in the centrifuge tube should dissolve during this period. After shaking, allow the layer of alkyl halides to separate (10–15 minutes at most). A fairly distinct top layer containing the products should have formed by this time.

CAUTION

***tert*-Butyl halides are volatile and should not be left in an open container any longer than necessary.**

Slowly remove most of the bottom aqueous layer with a Pasteur pipette and transfer it to a beaker. After waiting 10–15 seconds, remove the remaining lower layer in the centrifuge tube, including a small amount of the upper organic layer, to be certain that the organic layer is not contaminated by any water.

NOTE: For the following purification sequence, be certain to be well prepared. If you find that you are taking longer than 5 minutes to complete the entire sequence, you probably have affected your results adversely!

Using a dry Pasteur pipette, transfer the remainder of the alkyl halide layer into a small test tube (10 × 75 mm) containing about 0.05 g of solid sodium bicarbonate. As soon as the bubbling stops and a clear liquid is obtained, transfer it with a Pasteur pipette into a small, dry leak-proof vial, taking care not to transfer any solid.³ *Be sure the cap is screwed on tightly.* Do not store the liquid in a container with a cork or a rubber stopper, because these will absorb the halides. This sample can now be analyzed by as many of the methods in Experiment 22C as your instructor indicates. If possible, analyze the sample on the same day. When you have finished this procedure, return to Experiment 22A.

22C EXPERIMENT 22C*Analysis***PROCEDURE**

The ratio of 1-chlorobutane to 1-bromobutane, 2-chlorobutane to 2-bromobutane, or *tert*-butyl chloride to *tert*-butyl bromide must be determined. At your instructor's option, you may do this by one of three methods: gas chromatography, refractive index, or NMR spectroscopy. The products obtained from the reactions of 1-butanol and 2-butanol, however, cannot be analyzed by the refractive index

³ See footnote 2.

method (they contain pentane). The products obtained from the reaction of *tert*-butyl alcohol may be difficult to analyze by gas chromatography because the *tert*-butyl halides sometimes undergo elimination in the gas chromatograph.⁴

Gas Chromatography⁵ The instructor or a laboratory assistant may either make the sample injections or allow you to make them. In the latter case, your instructor will give you adequate instruction beforehand. A reasonable sample size is 2.5 μL . Inject the sample into the gas chromatograph and record the gas chromatogram. Determine the identity of the two peaks by comparing the retention times to the retention times of the pure compounds, which your instructor will provide.

Once the gas chromatogram has been obtained, determine the relative areas of the two peaks (Technique 22, Section 22.12). Remember to ignore the pentane peak. If the gas chromatograph has an integrator, it will report the areas. Triangulation is the preferred method of determining areas, if an integrator is not available. Record the percentages of alkyl chloride and alkyl bromide in the reaction mixture.

Nuclear Magnetic Resonance Spectroscopy

The instructor or a laboratory assistant will record the NMR spectrum of the reaction mixture.⁶ Submit a sample vial containing the mixture for this spectral determination. The spectrum will also contain integration of the important peaks (Technique 21, Nuclear Magnetic Resonance Spectroscopy).

If the substrate alcohol was 1-butanol, the resulting halide and pentane mixture will give rise to a complicated spectrum. Each alkyl halide will show a down-field triplet caused by the CH_2 group nearest the halogen. This triplet will appear farther downfield for the alkyl chloride than for the alkyl bromide. In a 60-MHz spectrum, these triplets will overlap, but one branch of each triplet will be available for comparison. Compare the integral of the *downfield* branch of the triplet for 1-chlorobutane with the *upfield* branch of the triplet for 1-bromobutane. The upper spectrum on the previous page provides an example. The relative heights of these integrals correspond to the relative amounts of each halide in the mixture.

If the substrate alcohol was 2-methyl-2-propanol, the resulting halide mixture will show two peaks in the NMR spectrum. Each halide will show a singlet because all the CH_3 groups are equivalent and are not coupled. In the reaction mixture, the upfield peak is due to *tert*-butyl chloride, and the downfield peak is caused by *tert*-butyl bromide. Compare the integrals of these peaks. The lower spectrum on the previous page provides an example. The relative heights of these integrals correspond to the relative amounts of each halide in the mixture.


REPORT

Record the percentages of alkyl chloride and alkyl bromide in the reaction mixture for each of the three alcohols. You need to share your data from the reaction with


⁴ *Note to the Instructor:* If pure samples of each product are available, check the assumption used here that the gas chromatograph responds equally to each substance. Response factors (relative sensitivities) are easily determined by injecting an equimolar mixture of products and comparing the peak areas.

⁵ *Note to the Instructor:* To obtain reasonable results for the gas chromatographic analysis of the *tert*-butyl halides, it may be necessary to supply the students with response factor correction (Technique 22, Section 22.13B).

⁶ It is difficult to determine the ratio of 2-chlorobutane to 2-bromobutane using nuclear magnetic resonance. This method requires at least a 90-MHz instrument. At 300 MHz, all downfield peaks are fully resolved.

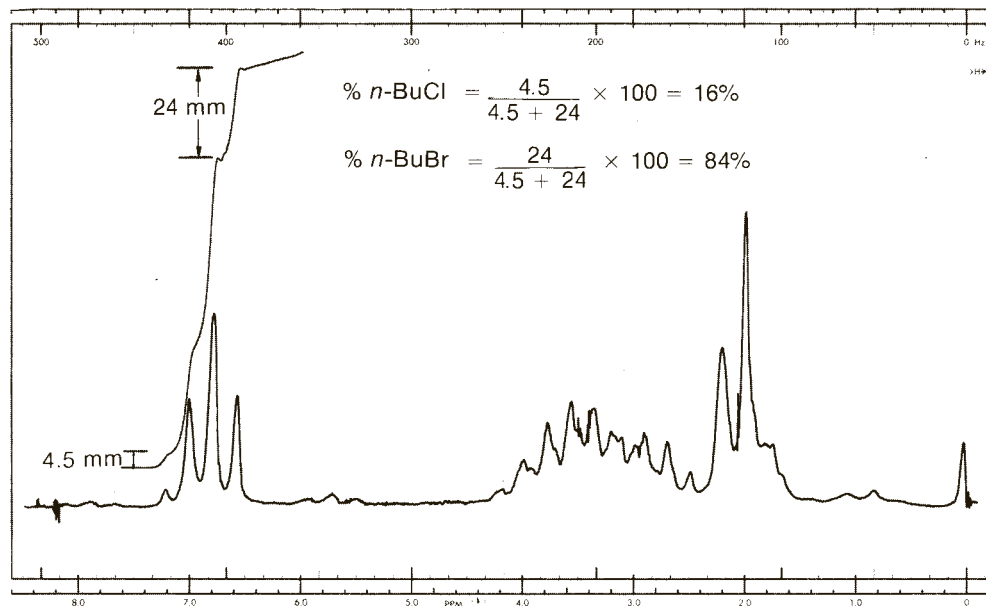

A 60-MHz NMR spectrum of 1-chlorobutane and 1-bromobutane, sweep width 250 Hz (no pentane in sample).

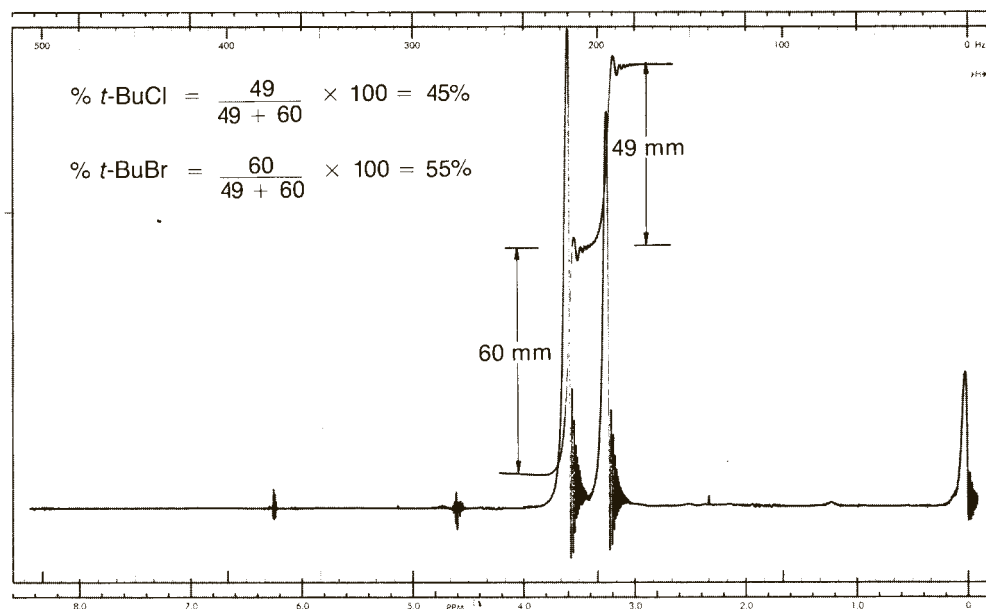

A 60-MHz spectrum of tert-butyl chloride and tert-butyl bromide, sweep width 250 Hz.

1-butanol or 2-butanol with other students in order to do this. The report must include the percentages of each alkyl halide determined by each method used in this experiment for the two alcohols you studied. On the basis of product distribution, develop an argument for which mechanism (S_N1 or S_N2) predominated for each of the three alcohols studied. The report should also include a discussion of which is the better nucleophile, chloride ion or bromide ion, based on the experimental results. All gas chromatograms and spectra should be attached to the report.

QUESTIONS

1. Draw complete mechanisms that explain the resultant product distributions observed for the reactions of *tert*-butyl alcohol and 1-butanol under the reaction conditions of this experiment.
2. Which is the better nucleophile in a protic solvent, chloride ion or bromide ion? Try to explain this in terms of the nature of the chloride ion and the bromide ion.
3. What is the principal organic by-product for each of these reactions?
4. A student left some alkyl halides (RCl and RBr) in an open container for several minutes. What happened to the composition of the halide mixture during that time? Assume that some liquid remains in the container.
5. What would happen if some of the solids in the nucleophile medium were not dissolved? How might this affect the outcome of the experiment?
6. What might have been the product ratios observed in this experiment if an aprotic solvent such as dimethyl sulfoxide had been used instead of water?
7. Explain the order of elution you observed while performing the gas chromatography for this experiment. What property of the product molecules seems to be the most important in determining relative retention times?
8. When you calculate the percentage composition of the product mixture, exactly what kind of "percentage" (i.e., volume percent, weight percent, mole percent) are you dealing with?
9. When a pure sample of *tert*-butyl bromide is analyzed by gas chromatography, two components are usually observed. One of them is *tert*-butyl bromide and the other one is a decomposition product. As the temperature of the injector is increased, the amount of the decomposition product increases and the amount of *tert*-butyl bromide decreases.
 - a. What is the structure of the decomposition product?
 - b. Why does the amount of decomposition increase with increasing temperature?
 - c. Why does *tert*-butyl bromide decompose much more easily than *tert*-butyl chloride?


Synthesis of *n*-Butyl Bromide and *t*-Pentyl Chloride

Synthesis of alkyl halides


Extraction

Simple distillation

The synthesis of two alkyl halides from alcohols is the basis for these experiments. In the first experiment, a primary alkyl halide *n*-butyl bromide is prepared as shown in Equation 1.

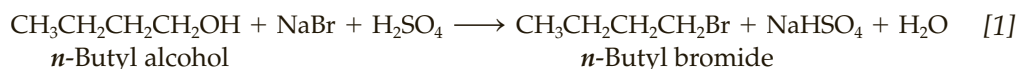


In the second experiment, a tertiary alkyl halide *t*-pentyl chloride (*t*-amyl chloride) is prepared as shown in Equation 2.

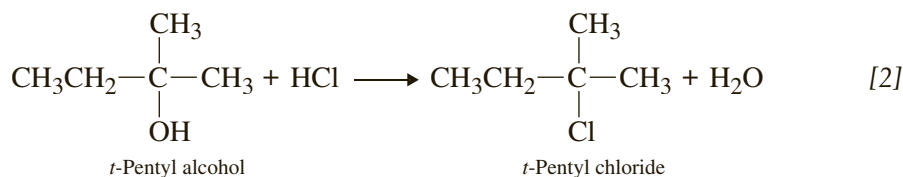


These reactions provide an interesting contrast in mechanisms. The *n*-butyl bromide synthesis proceeds by an $\text{S}_{\text{N}}2$ mechanism, whereas *t*-pentyl chloride is prepared by an $\text{S}_{\text{N}}1$ reaction.


n-Butyl Bromide

The primary alkyl halide *n*-butyl bromide can be prepared easily by allowing *n*-butyl alcohol to react with sodium bromide and sulfuric acid by Equation 1. The sodium bromide reacts with sulfuric acid to produce hydrobromic acid.

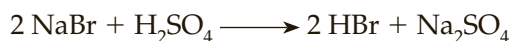

Excess sulfuric acid serves to shift the equilibrium and thus to speed the reaction by producing a higher concentration of hydrobromic acid. The sulfuric acid also protonates the hydroxyl group of *n*-butyl alcohol so that water is displaced. The acid also protonates the water as it is produced in the reaction and deactivates it as a nucleophile. Deactivation of water keeps the alkyl halide from being converted back to the alcohol by nucleophilic attack of water. The reaction of the primary substrate proceeds via an $\text{S}_{\text{N}}2$ mechanism.

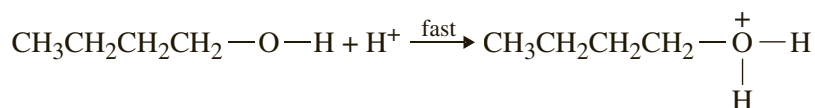


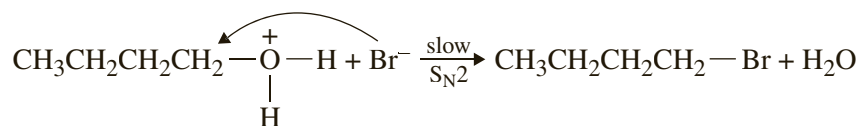

During the isolation of the *n*-butyl bromide, the crude product is washed with sulfuric acid, water, and sodium bicarbonate to remove any remaining acid or *n*-butyl alcohol.

t-Pentyl Chloride

The tertiary alkyl halide can be prepared by allowing *t*-pentyl alcohol to react with concentrated hydrochloric acid according to Equation 2. The reaction is accomplished simply by shaking the two reagents in a sealed conical vial. As the reaction proceeds, the insoluble alkyl halide product forms on upper phase. The reaction of the tertiary substrate occurs via an $\text{S}_{\text{N}}1$ mechanism.

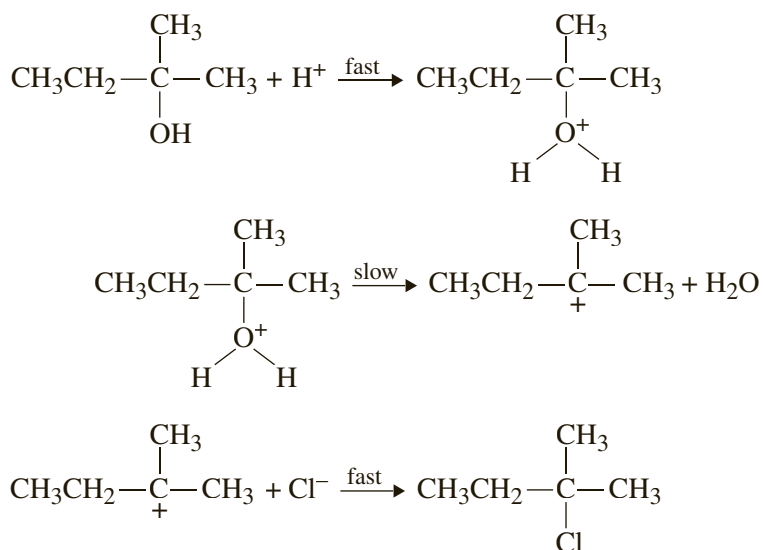

A small amount of alkene, 2-methyl-2-butene, is produced as a by-product in this reaction. If sulfuric acid had been used as it was for *n*-butyl bromide, a much larger amount of this alkene would have been produced.

REQUIRED READING

Review: Techniques 1 through 7
Techniques 12, 13, and 14

SPECIAL INSTRUCTIONS

CAUTION

Take special care with concentrated sulfuric acid: It causes severe burns.



As your instructor indicates, perform either the *n*-butyl bromide or the *t*-pentyl chloride procedure or both.

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous solutions produced in this experiment in the container for aqueous waste.

NOTE TO THE INSTRUCTOR

The temperature during the distillations in the various parts of Experiment 23 may be monitored with either a thermometer or a stainless steel temperature probe. If a stainless steel probe is used, it must be used in conjunction with either a digital thermometer or one of the Vernier devices (see Technique 13, Section 13.4, and Technique 14, Figure 14.12).

For the boiling point determination, we prefer the Semimicroscale Direct Method described in Technique 13, Section 13.2. The best way to perform this method is to use a digital thermometer with a stainless steel probe (see Technique 13, Section 13.4 and Figure 13.7).

23A EXPERIMENT 23A

n-Butyl Bromide

PROCEDURE

Preparation of n-Butyl Bromide


Using an automatic pipette or a dispensing pump, place 1.4 mL of *n*-butyl alcohol (1-butanol, MW = 74.1) in a preweighed 10-mL round-bottom flask. Reweigh the flask to determine the exact weight of the alcohol. Add 2.4 g of sodium bromide and 2.4 mL of water. Cool the mixture in an ice bath and slowly add 2.0 mL of concentrated sulfuric acid dropwise using a Pasteur pipette. Add a magnetic stirring bar and assemble the reflux apparatus and trap shown in the figure. The trap absorbs the hydrogen bromide gas evolved during the reaction period. While stirring, heat the mixture to a gentle boil (aluminum block temperature about 145°C) for 60–75 minutes.

Extraction

Remove the heat source and allow the apparatus to cool until you can disconnect the round-bottom flask without burning your fingers.

NOTE: Do not let the reaction mixture cool to room temperature. Complete the operations in this paragraph as quickly as possible. Otherwise, salts may precipitate, making this procedure more difficult to perform.

The *n*-butyl bromide layer should be on top. If the reaction is not yet complete, the remaining *n*-butyl alcohol will sometimes form a *second organic layer* on top of the *n*-butyl bromide layer. Treat both organic layers as if they were one. Remove and discard as much of the aqueous (bottom) layer as possible using a Pasteur pipette, but do not remove any of the organic layer (or layers). Ignore the salts during this

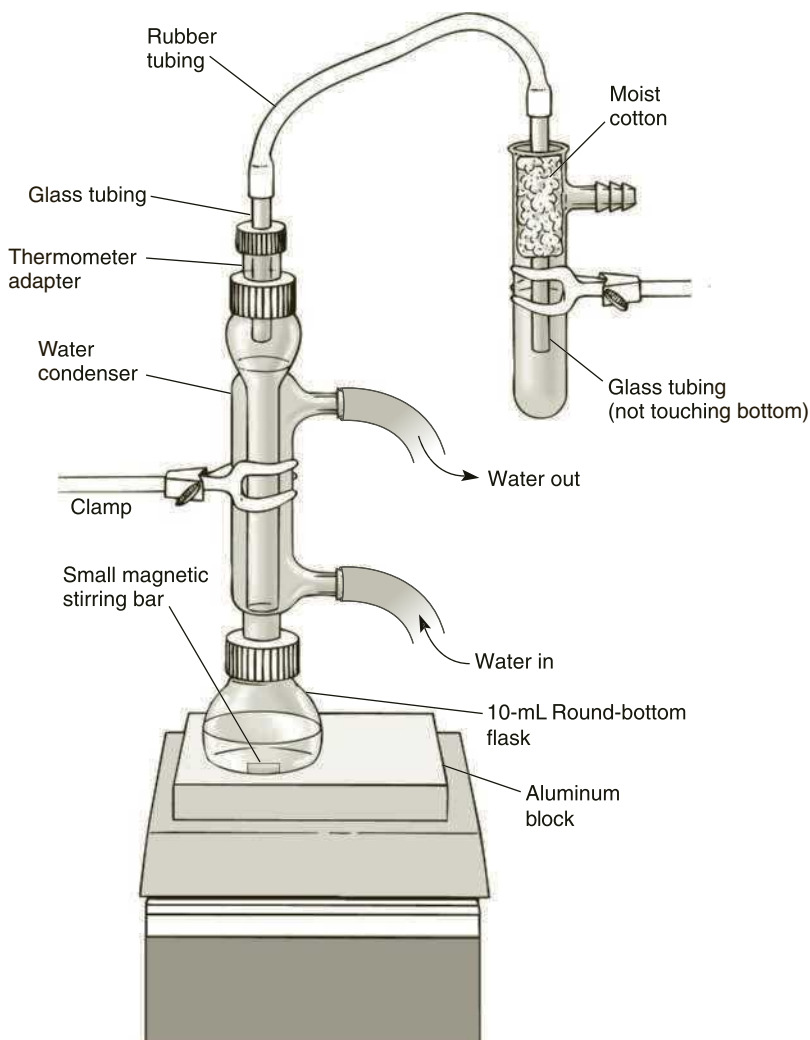

Apparatus for Experiment 23A, *n*-butyl bromide.

separation. If they are drawn into the pipette, treat them as part of the aqueous layer. Transfer the remaining liquid to a 5-mL conical vial. Remove and discard any aqueous layer remaining in the conical vial.

The organic and aqueous layers should separate as described in the following instructions. However, to make sure that you do not discard the wrong layer, it would be a good idea to add a drop of water to any aqueous layer you plan to discard. If a drop of water dissolves in the liquid, you can be confident that it is an aqueous layer. Add 2 mL of 9M H_2SO_4 to the conical vial. Cap the vial and shake it *gently*, venting occasionally. Allow the layers to separate. Because any remaining *n*-butyl alcohol is extracted by the H_2SO_4 solution, there should now be only one organic layer. The organic layer should be the top layer. Remove and discard the aqueous (bottom) layer.

Add 2 mL of water to the vial. Cap the vial and shake it *gently*, venting occasionally. Allow the layers to separate. This time, the organic layer should be the bottom layer. The bottom layer may form into a globule (ball) instead of separating cleanly. Use a microspatula to prod the ball gently into the bottom of the vial. Using a Pasteur pipette, transfer the bottom layer (or globule) to a clean 5-mL conical vial. Add 2 mL of saturated aqueous sodium bicarbonate solution, a little at a time, while stirring. Cap the vial and shake it vigorously for 1 minute, venting frequently to relieve any pressure that is produced. Allow the layers to separate and then carefully transfer the

lower alkyl halide layer to a dry 3-mL conical vial using a dry Pasteur pipette. Dry the liquid over granular anhydrous sodium sulfate (Technique 12, Section 12.9).

Distillation

When the solution is dry, transfer it to a clean, dry, 3-mL vial using a Pasteur pipette and distill it (aluminum block about 140°C) using a clean, dry Hickman still (Technique 14, Figure 14.5). Each time the Hickman head becomes full, transfer the distillate to a preweighed conical vial using a Pasteur pipette.


When the distillation is complete (one or two drops remaining), weigh the vial, calculate the percentage yield, and determine a microscale boiling point (Technique 13, Section 13.2). Determine the infrared spectrum of the product using salt plates (Technique 25, Section 25.2). Submit the remainder of your sample in a properly labeled vial, along with the infrared spectrum, when you submit your report to the instructor.

23B EXPERIMENT 23B

n-Butyl Bromide (Semimicroscale Procedure)

PROCEDURE

Follow the procedure given in Experiment 23A, except double the amounts of *all* reagents. Use a 25-mL round-bottom flask rather than a 10-mL round-bottom flask for running the reaction. For the separation and extraction procedures, use a screw-cap centrifuge tube in place of a 5-mL conical vial. Distill the crude *n*-butyl bromide using an apparatus similar to the semimicroscale apparatus for a simple distillation (see Figure 14.10). Make the following changes: Use a 5-mL conical vial as the distilling flask and collect the distillate in a preweighed 3-mL conical vial rather than a graduated cylinder. The bulb of the thermometer or the bottom of the temperature probe must be placed below the side arm or it will not read the temperature correctly. All the glassware must be dry. Use a boiling stone or magnetic spin vane to prevent bumping. Collect the material that boils between 85°C and 102°C.

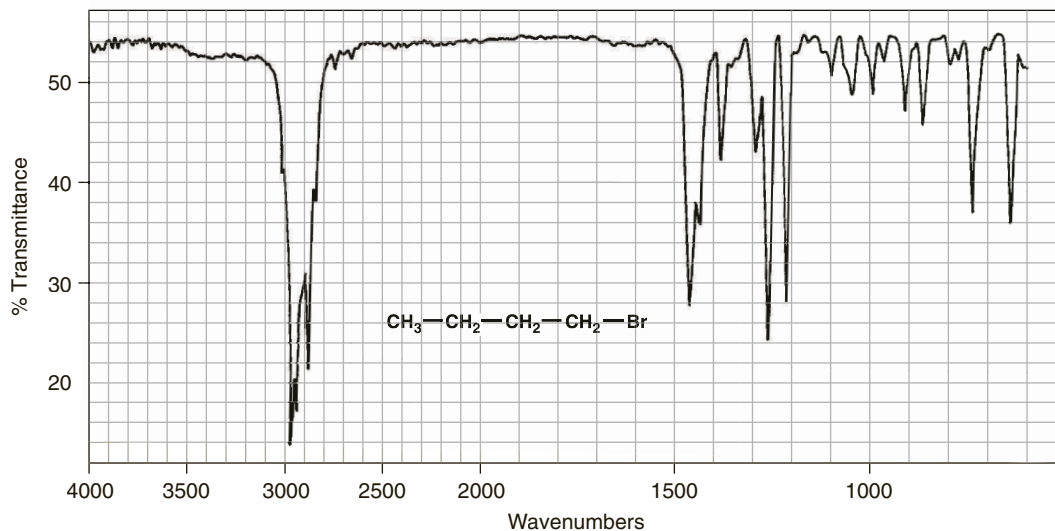

Infrared spectrum of *n*-butyl bromide (neat).

23C EXPERIMENT 23C*t*-Pentyl Chloride (Microscale Procedure)**PROCEDURE**

NOTE: In the following procedures, it may be difficult to see the interfaces between layers because the refractive index of the product will be similar to the refractive indices of the extraction solvents.

Preparation of *t*-Pentyl Chloride

Using an automatic pipette or a dispensing pump, place 1.0 mL of *t*-pentyl alcohol (2-methyl-2-butanol, MW = 88.2) in a preweighed 5-mL conical vial. Reweigh the vial to determine the exact weight of alcohol delivered.

NOTE: Before shaking the conical vial vigorously in the next step, be sure that the capped vial does not leak. If it does leak, use a Pasteur pipette to mix the two layers. Draw up as much liquid as possible into the Pasteur pipette and then expel the liquid rapidly back into the conical vial. Continue this mixing for 3–4 minutes.

Add 2.5 mL of concentrated hydrochloric acid, cap the vial, and shake it vigorously for 1 minute. After shaking the vial, loosen the cap and vent the vial. Recap the vial and shake it for 3 minutes more, venting occasionally. Allow the mixture to stand in the vial until the layer of alkyl halide product separates. The *t*-pentyl chloride ($d = 0.865$ g/mL) should be the top layer, but be sure to verify this carefully by observation as you add a few drops of hydrochloric acid.

With a Pasteur pipette, separate the layers by placing the tip of the pipette squarely into the bottom of the vial and removing the lower (aqueous) layer. Discard the aqueous layer. (Are you sure which one it is?)

Extraction

Carry out the operations in this paragraph as rapidly as possible because the *t*-pentyl chloride is unstable in water and aqueous bicarbonate solution. It is easily hydrolyzed back to the alcohol. Be sure everything you need is at hand. In each of the following steps, the organic layer should be on top; however, you should add a few drops of water to make sure. Wash the organic layer by adding 1 mL of water to the conical vial. Shake the mixtures for a few seconds and then allow the layers to re-form. Once again, separate the layers using a Pasteur pipette and discard the aqueous layer after making certain that you saved the proper layer. Add a 1-mL portion of 5% aqueous sodium bicarbonate to the organic layer. *Gently* mix the two phases in the vial with a stirring rod until they are thoroughly mixed. Now cap the vial and shake it gently for 1 minute, venting occasionally. After this, vigorously shake the vial for another 30 seconds, venting occasionally. Discard the aqueous layer and transfer the organic layer to a dry conical vial with a dry Pasteur pipette.

Dry the crude *t*-pentyl chloride over granular anhydrous sodium sulfate (Technique 12, Section 12.9).

Distillation

When the solution is dry (it should be clear), carefully separate the alkyl halide from the drying agent with a Pasteur pipette and transfer it to a clean, dry 3-mL conical vial. Add a microporous boiling stone and distill the crude *t*-pentyl chloride (Technique 14, Figure 14.5, or, if possible, Figure 14.7B).


Using a Pasteur pipette, transfer the product to a dry, preweighed conical vial, weigh it, and calculate the percentage yield. Determine a boiling point for the product using a microscale boiling-point determination (Technique 13, Section 13.2). Determine the infrared spectrum of the alkyl halide using salt plates (Technique 25, Section 25.2). Submit the remainder of your sample in a properly labeled vial, along with the infrared spectrum, when you submit your report to the instructor.

23D EXPERIMENT 23D

t-Pentyl Chloride (Semimicroscale Procedure)

PROCEDURE

Follow the procedure as written in Experiment 23C, except double the amounts of *all* reagents. Use a 15-mL screw-cap centrifuge tube instead of a 5-mL conical vial for running the reaction and performing the extractions. Distill the crude *t*-pentyl chloride using an apparatus similar to the semi-microscale apparatus for a simple distillation (see Technique 14, Figure 14.10). Make the following changes: Use a 5-mL conical vial as the distilling flask and collect the distillate in a preweighed 3-mL conical vial rather than a graduated cylinder. The bulb of the thermometer or the bottom of the temperature probe must be placed below the side arm or it will not read the temperature correctly. All the glassware must be dry. Use a boiling stone or magnetic spin vane to prevent bumping. Collect the material that boils between 80°C and 84°C.

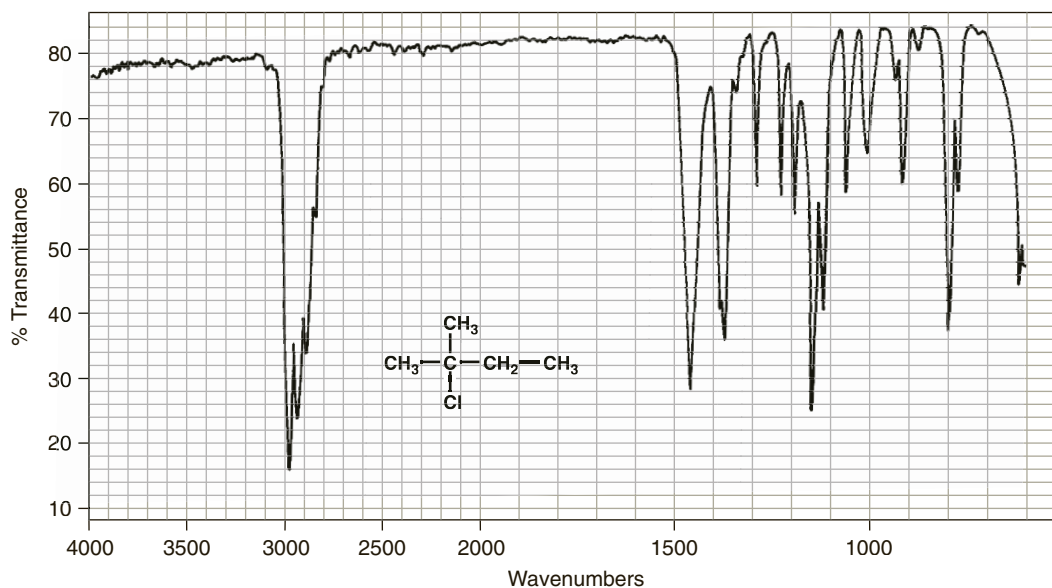

Infrared spectrum of *tert*-pentyl chloride (neat).

23E EXPERIMENT 23E*t*-Pentyl Chloride (Macroscale Procedure)**PROCEDURE***Preparation of t-Pentyl Chloride*

In a 125-mL separatory funnel, place 10.0 mL of *tert*-pentyl alcohol (2-methyl-2-butanol, $MW = 88.2$, $d = 0.805$ g/mL) and 25 mL of concentrated hydrochloric acid ($d = 1.18$ g/mL). Do not stopper the funnel. Gently swirl the mixture in the separatory funnel for about 1 minute. After this period of swirling, stopper the separatory funnel and carefully invert it. Without shaking the separatory funnel, immediately open the stopcock to release the pressure. Close the stopcock, shake the funnel several times, and again release the pressure through the stopcock (Technique 12, Section 12.7). Shake the funnel for 2–3 minutes, with occasional venting. Allow the mixture to stand in the separatory funnel until the two layers have completely separated. The *tert*-pentyl chloride ($d = 0.865$ g/mL) should be the top layer, but be sure to verify this by adding a few drops of water. The water should dissolve in the lower (aqueous) layer. Drain and discard the lower layer.

Extraction

The operations in this paragraph should be done as rapidly as possible because the *tert*-pentyl chloride is unstable in water and sodium bicarbonate solution. It is easily hydrolyzed back to the alcohol. In each of the following steps, the organic layer should be on top; however, you should add a few drops of water to make sure. Wash (swirl and shake) the organic layer with 10 mL of water. Separate the layers and discard the aqueous phase after making certain that the proper layer has been saved. Add a 10-mL portion of 5% aqueous sodium bicarbonate to the separatory funnel. Gently swirl the funnel (unstoppered) until the contents are thoroughly mixed. Stopper the funnel and carefully invert it. Release the excess pressure through the stopcock. Gently shake the separatory funnel, with frequent release of pressure. After this, vigorously shake the funnel, again with release of pressure, for about 1 minute. Allow the layers to separate and drain the lower aqueous layer. Wash (swirl and shake) the organic layer with one 10-mL portion of water and again drain the lower aqueous layer. Transfer the organic layer to a small, dry Erlenmeyer flask by pouring it from the top of the separatory funnel. Dry the crude *tert*-pentyl chloride over 1.0 g of anhydrous calcium chloride until it is clear (Technique 12, Section 12.9). Swirl the alkyl halide with the drying agent to aid the drying.

Distillation

Transfer the clear liquid to a dry 25-mL round-bottom flask using a Pasteur pipette. Add a boiling stone and distill the crude *tert*-pentyl chloride in a dry apparatus (Technique 14, Section 14.4, Figure 14.11). Collect the pure *tert*-pentyl chloride in a receiver cooled in ice. Collect the material that boils between 78°C and 84°C. Weigh the product, calculate the percentage yield, and determine the boiling point (Technique 13, Section 13.2). Determine the infrared spectrum of the product using salt plates (Technique 25, Section 25.2). Submit the remainder of your sample in a properly labeled vial, along with the infrared spectrum, when you submit your report to the instructor.

QUESTIONS

n-Butyl Bromide

1. What are the formulas of the salts that precipitate when the reaction mixture is cooled?
2. Why does the alkyl halide layer switch from the top layer to the bottom layer at the point where water is used to extract the organic layer?
3. An ether and an alkene are formed as by-products in this reaction. Draw the structures of these by-products and give mechanisms for their formation.
4. Aqueous sodium bicarbonate was used to wash the crude *n*-butyl bromide.
 - a. What was the purpose of this wash? Give equations.
 - b. Why would it be undesirable to wash the crude halide with aqueous sodium hydroxide?
5. Look up the density of *n*-butyl chloride (1-chlorobutane). Assume that this alkyl halide was prepared instead of the bromide. Decide whether the alkyl chloride would appear as the upper or the lower phase at each stage of the separation procedure: after the reflux, after the addition of water, and after the addition of sodium bicarbonate.
6. Why must the alkyl halide product be dried carefully with anhydrous sodium sulfate before the distillation? (*Hint*: See Technique 15, Section 15.7.)

t-Pentyl Chloride

1. Aqueous sodium bicarbonate was used to wash the crude *t*-pentyl chloride.
 - a. What was the purpose of this wash? Give equations.
 - b. Why would it be undesirable to wash the crude halide with aqueous sodium hydroxide?
2. Some 2-methyl-2-butene may be produced in the reaction as a by-product. Give a mechanism for its production.
3. How is unreacted *t*-pentyl alcohol removed in this experiment? Look up the solubility of the alcohol and the alkyl halide in water.
4. Why must the alkyl halide product be dried carefully with anhydrous sodium sulfate before the distillation? (*Hint*: See Technique 15, Section 15.7.)
5. Will *t*-pentyl chloride (2-chloro-2-methylbutane) float on the surface of water? Look up its density in a handbook.


4-Methylcyclohexene

Preparation of an alkene

Dehydration of an alcohol

Distillation

Bromine and permanganate tests for unsaturation

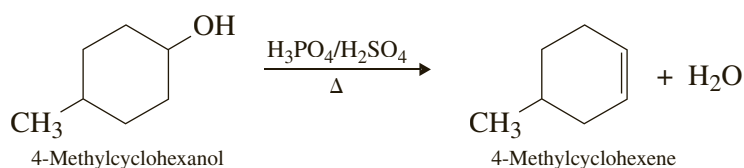


Alcohol dehydration is an acid-catalyzed reaction performed by strong, concentrated mineral acids such as sulfuric and phosphoric acids. The acid protonates the alcoholic hydroxyl group, permitting it to dissociate as water. Loss of a proton from the intermediate (elimination) brings about an alkene. Because sulfuric acid often causes extensive charring in this reaction, phosphoric acid, which is comparatively free of this problem, is a better choice. To make the reaction proceed faster, however, a minimal amount of sulfuric acid will also be used.

The equilibrium that attends this reaction will be shifted in favor of the product by distilling it from the reaction mixture as it is formed. The 4-methylcyclohexene (bp 101–102°C) will codistill with the water that is also formed. By continuously removing the products, one can obtain a high yield of 4-methylcyclohexene. Because the starting material, 4-methylcyclohexanol, also has a somewhat low boiling point (bp 171–173°C), the distillation must be done carefully so that the alcohol does not also distill.

Unavoidably, a small amount of phosphoric acid codistills with the product. It is removed by washing the distillate mixture with a saturated sodium chloride solution. This step also partially removes the water from the 4-methylcyclohexene layer; the drying process will be completed by allowing the product to stand over anhydrous sodium sulfate.

Compounds containing double bonds react with a bromine solution (red) to decolorize it. Similarly, they react with a solution of potassium permanganate (purple) to discharge its color and produce a brown precipitate (MnO₂). These reactions are often used as qualitative tests to determine the presence of a double bond in an organic molecule (see Experiment 57). Both tests will be performed on the 4-methylcyclohexene formed in this experiment.

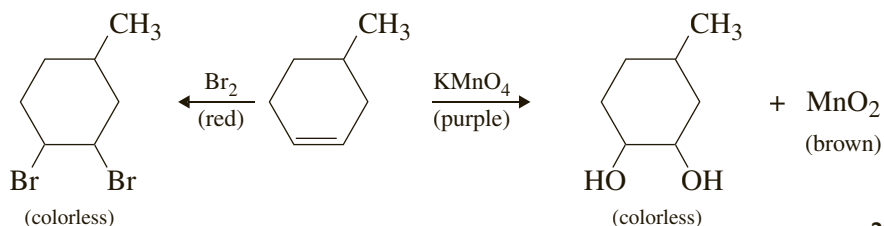

REQUIRED READING

Review: Techniques 1, 2, 3, 5, and 6

Technique 12 Extractions, Separations, and Drying Agents,
Sections 12.5 and 12.9

New: Technique 14 Simple Distillation

If performing the optional boiling point or infrared spectroscopy, also read:

Technique 13 Physical Constants of Liquids: The Boiling Point
and Density

Technique 25 Infrared Spectroscopy

SPECIAL INSTRUCTIONS

Phosphoric and sulfuric acids are corrosive. Do not allow either acid to touch your skin.

If you must store the 4-methylcyclohexene, it is essential that the sample be stored in leak-proof vials. We have found that GC-MS vials work much better for this purpose than a conical vial with a septum liner and cap.

SUGGESTED WASTE DISPOSAL

Any organic residues should be discarded in an organic waste container designated for the disposal of *nonhalogenated* wastes. Discard the solutions that remain after the bromine test for unsaturation in an organic waste container designated for the disposal of *halogenated* wastes. The solutions that remain after the potassium permanganate test should be discarded into a waste container specifically marked for the disposal of heavy-metal wastes. Aqueous solutions should be placed in the container designated for that purpose.

NOTE TO THE INSTRUCTOR

Amberlyst-15 ion exchange resin (sulfonic acid groups) may be used in place of the phosphoric and sulfuric acids.¹ There is less charring with the resin. Use 0.2 g of the resin and heat more slowly, increasing the reaction time to about 45 minutes. When measuring the resin (little balls), use a measuring spoon with a depression; the spheres roll off a flat spatula, and static charges sometimes complicate the weighing problem. Provide a waste container for the spent resin.

The temperature during the distillation in Experiment 24B may be monitored with either a thermometer or a stainless steel temperature probe. If a stainless steel probe is used, it must be used in conjunction with either a digital thermometer or one of the Vernier devices (see Technique 13, Section 13.4, and Technique 14, Figure 14.12).

¹ Moeur, H. P.; Swatik, S. A.; and Pinnell, R. P. Microscale Dehydration of Cyclohexanol Using a Macroreticular Cation Exchange Resin. *J. Chem. Edu.* **1997**, *74* (Jul); 833.

For the boiling point determination, we prefer the Semimicroscale Direct Method described in Technique 13, Section 13.2. The best way to perform this method is to use a digital thermometer with a stainless steel probe (see Technique 13, Section 13.4 and Figure 13.7).

24A EXPERIMENT 24A

4-Methylcyclohexene (Microscale Procedure)

PROCEDURE

Apparatus Assembly

Place 1.5 mL of 4-methylcyclohexanol ($MW = 114.2$) in a tared 5-mL conical vial and reweigh the vial to determine an accurate weight for the alcohol. Add 0.40 mL of 85% phosphoric acid and six drops of concentrated sulfuric acid to the vial. Mix the liquids thoroughly using a glass stirring rod and add a boiling stone. Assemble a distillation apparatus as shown in Figure 14.5 and use a water-cooled condenser. It is recommended that you include the drying tube, filled with calcium chloride, as an odor-control measure.

Dehydration

Start circulating the cooling water in the condenser and heat the mixture until the product begins to distill (aluminum block or sand bath set to about 160–180°C). If you use an aluminum block for heating, place aluminum collars around the conical vial. The heating should be regulated so that the distillation requires about 30–45 minutes, heating slowly at the beginning.

During the distillation, use a Pasteur pipette to remove the distillate from the well of the Hickman head when necessary. You must remove the condenser when performing this experiment, unless you have a Hickman head with a side port. In that case, you can remove the distillate through the side port without removing the condenser. Transfer the distillate to a clean, dry, 3-mL conical vial, which should be capped except when liquid is being added. Continue the distillation until no more liquid is being distilled. This can be best determined by observing when boiling in the conical vial has ceased. Also, the volume of liquid remaining in the vial should be about 0.5 mL when distillation is complete.

When distillation is complete, remove as much distillate as possible from the Hickman head and transfer it to the 3-mL conical vial. Then, using a Pasteur pipette with the tip slightly bent, rinse the sides of the inside wall of the Hickman head with 1.0 mL of saturated sodium chloride. Do this thoroughly so that as much liquid as possible is washed down into the well of the Hickman head. Transfer this liquid to the 3-mL conical vial.

Isolation and Drying of Product

Allow the layers to separate and remove the bottom aqueous layer. Using a dry Pasteur pipette, transfer the organic layer to a small test tube, and dry it over granular anhydrous sodium sulfate (Technique 12, Section 12.9). Place a stopper in the test tube and set it aside for 10–15 minutes to remove the last traces of water. Carefully

transfer as much distillate as possible to a small tared conical vial with a cap or a GC-MS vial. Weigh the product ($MW = 96.2$) and calculate the percentage yield.

Boiling-Point Determination and Spectroscopy²

At the instructor's option, determine an accurate boiling point on your sample using the microboiling-point method (Technique 13, Section 13.2), and obtain the infrared spectrum of 4-methylcyclohexene (Technique 25, Section 25.2 or Section 25.3). Because 4-methylcyclohexene is so volatile, you must work quickly to obtain a good spectrum using sodium chloride plates. A better method is to use silver chloride plates with the orientation shown in Figure 25.3B. Compare the spectrum with the one shown in this experiment. After performing the tests at the end of this experiment, submit your sample, along with the report, to the instructor.

24B EXPERIMENT 24B

4-Methylcyclohexene (Semimicroscale Procedure)

PROCEDURE

Apparatus Assembly

Assemble a distillation apparatus as shown in Figure 14.10, but insert a water condenser as shown in Figure 14.11 and insert the thermometer or temperature probe a bit lower than shown so that the mercury bulb or end of the temperature probe is inside the ground-glass joint. Place 4.0 mL of 4-methylcyclohexanol ($MW = 114.2$) in a 10-mL graduated cylinder and weigh it. Transfer the alcohol to the 10-mL round-bottom flask and reweigh the graduated cylinder to determine an accurate weight for the alcohol. Add 1.0 mL of 85% phosphoric acid and 16 drops of concentrated sulfuric acid to the alcohol already in the round-bottom flask. Without delay, mix the contents of the round-bottom flask thoroughly by swirling the liquids in the flask. Add a corundum (black) boiling stone and reconnect the flask to the distillation apparatus.

Dehydration

Start circulating cooling water in the condenser and heat the mixture until the product begins to distill (aluminum block approximately 180°C). The heating should be regulated so that the distillation requires about 45 minutes to be completed. Distill slowly and monitor (record) the temperature while the liquid distills; the products should distill over a range of about $90\text{--}105^{\circ}\text{C}$. Continue distillation until no more liquid is collected. This can best be determined by observing when the temperature drops. You should collect about 3 mL of distillate. Approximately 1 mL of a dark brown liquid will remain undistilled in the round-bottom flask when the distillation is complete.

² An option is to submit the student samples to analysis by gas chromatography–mass spectrometry. If GC-MS is used, students can observe that several products besides 4-methylcyclohexene are formed in the reaction. Discussion of the mechanism of rearrangement should then be included.


When the distillation is complete, open the joint between the distillation head and the condenser, and then tilt the condenser so that any liquid trapped in the joints of the condenser will drip into the collection vial. If droplets of distillate are evident clinging to the inside tube of the condenser, wash them down with a small amount of saturated aqueous sodium chloride solution.

Isolation and Drying of the Product

Using a Pasteur pipette, transfer the distillate from the conical vial to a 15-mL capped centrifuge tube, and add approximately 2.0 mL of a saturated sodium chloride solution. Cap the tube and gently invert it several times with venting. Allow the layers to separate and remove the lower aqueous layer with a Pasteur pipette. Transfer the organic layer to a small test tube or Erlenmeyer flask, and dry it over granular anhydrous sodium sulfate (Technique 12, Section 12.9). Place a stopper in the test tube (or Erlenmeyer flask), and set it aside for 10–15 minutes to remove the last traces of water. Transfer the dry distillate to a small tared conical vial with a cap or a GC-MS vial. Weigh the product ($MW = 96.2$) and calculate the percentage yield.


Boiling-Point Determination and Spectroscopy³

At the instructor's option, determine a more accurate boiling point on your sample using the micro boiling-point method (Technique 13, Section 13.2) and obtain the infrared spectrum of 4-methylcyclohexene (Technique 25, Section 25.2 or Section 25.3). Because 4-methylcyclohexene is so volatile, you must work quickly to obtain a good spectrum using sodium chloride plates. A better method is to use silver chloride plates with the orientation shown in Figure 25.3B. Compare the spectrum with the one shown in this experiment. After performing the following tests, submit your sample, along with the report, to the instructor.

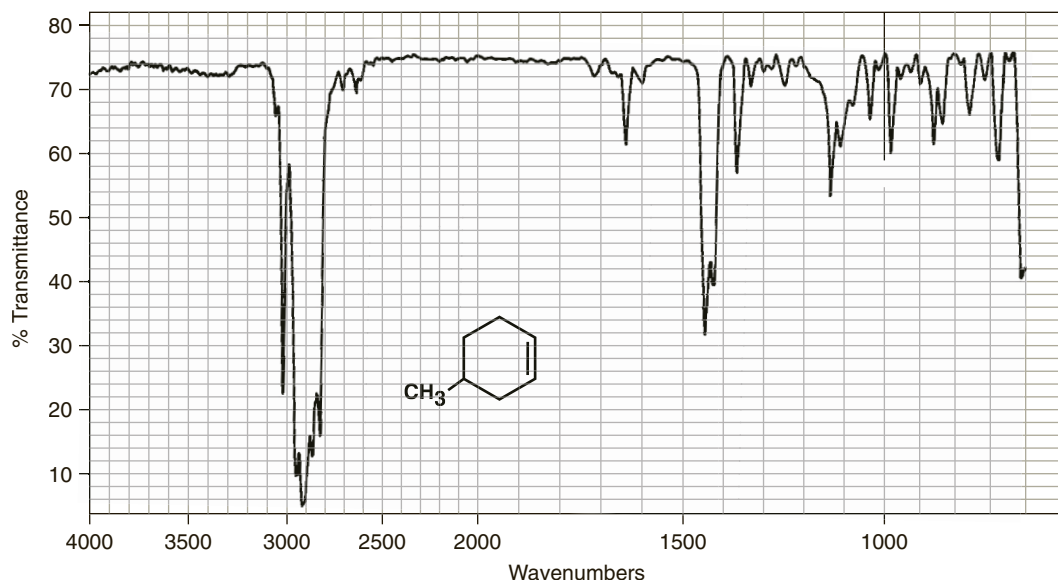

Infrared spectrum of 4-methylcyclohexene (neat).

³ See footnote 2 in Experiment 24A.

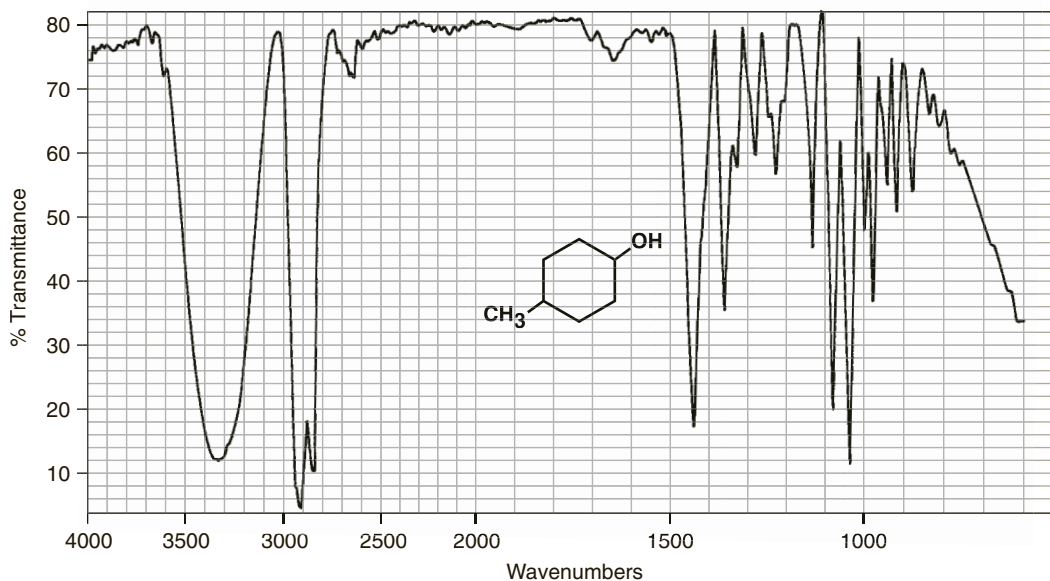

Infrared spectrum of 4-methylcyclohexanol (neat).

UNSATURATION TESTS

Place four to five drops of 4-methylcyclohexanol in each of two small test tubes. In each of another pair of small test tubes, place four to five drops of the 4-methylcyclohexene you prepared. Do not confuse the test tubes. Take one test tube from each group, and add a solution of bromine in carbon tetrachloride or methylene chloride, drop by drop, to the contents of the test tube, until the red color is no longer discharged or until you have added 20 drops. Record the result in each case, including the number of drops required. Test the remaining two test tubes in a similar fashion with a solution of potassium permanganate. Because aqueous potassium permanganate is not miscible with organic compounds, you will have to add about 0.3 mL of 1,2-dimethoxyethane to each test tube before making the test. Record your results and explain them.


QUESTIONS

- Outline a mechanism for the dehydration of 4-methylcyclohexanol catalyzed by phosphoric acid.
- What major alkene product is produced by dehydrating the following alcohols?
 - Cyclohexanol
 - 1-Methylcyclohexanol
 - 2-Methylcyclohexanol
 - 2,2-Dimethylcyclohexanol
 - 1,2-Cyclohexanediol (*Hint*: Consider keto–enol tautomerism.)
- Compare and interpret the infrared spectra of 4-methylcyclohexene and 4-methylcyclohexanol.
- Identify the C-H out-of-plane bending vibrations in the infrared spectrum of 4-methylcyclohexene. What structural information can be obtained from these bands?
- In this experiment, 1.0 mL of saturated sodium chloride is used to rinse the Hickman head after the initial distillation. Why is saturated sodium chloride, rather than pure water, used for this procedure and the subsequent washing of the organic layer?

Fats and Oils

In the normal human diet, about 25% to 50% of the caloric intake consists of fats and oils. These substances are the most concentrated form of food energy in our diet. When metabolized, fats produce about 9.5 kcal of energy per gram. Carbohydrates and proteins produce less than half this amount. For this reason, animals tend to build up fat deposits as a reserve source of energy. They do this, of course, only when their food intake exceeds their energy requirements. In times of starvation, the body metabolizes these stored fats. Even so, some fats are required by animals for bodily insulation and as a protective sheath around some vital organs.

The constitution of fats and oils was first investigated by the French chemist Chevreul from 1810 to 1820. He found that when fats and oils were hydrolyzed, they gave rise to several "fatty acids" and the trihydroxylic alcohol glycerol. Thus, fats and oils are **esters** of glycerol, called **glycerides** or **acylglycerols**. Because glycerol has three hydroxyl groups, it is possible to have mono-, di-, and triglycerides. Fats and oils are predominantly triglycerides (triacylglycerols), constituted as follows:

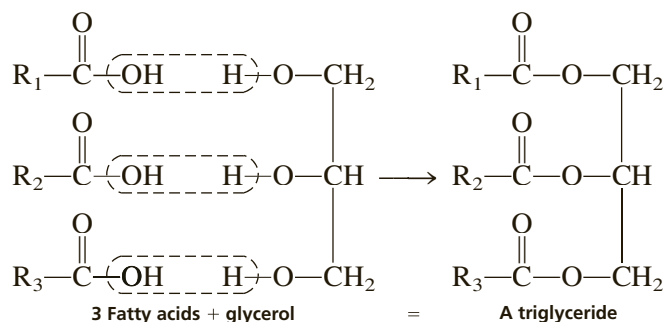

Thus, most fats and oils are esters of glycerol, and their differences result from the differences in the fatty acids with which glycerol may be combined. The most common fatty acids have 12, 14, 16, or 18 carbons, although acids with both lesser and greater numbers of carbons are found in several fats and oils. These common fatty acids are listed in Table 1 along with their structures. As you can see, these acids are both saturated and unsaturated. The saturated acids tend to be solids, whereas the unsaturated acids are usually liquids. This circumstance also extends to fats and oils. Fats are made up of fatty acids that are most saturated, whereas oils are primarily composed of fatty acid portions that have greater numbers of double bonds. In other words, unsaturation lowers the melting point. Fats (solids) are usually obtained from animal sources, whereas oils (liquids) are commonly obtained from vegetable sources. Therefore, vegetable oils usually have a higher degree of unsaturation.

About 20 to 30 fatty acids are found in fats and oils, and it is not uncommon for a given fat or oil to be composed of as many as 10 to 12 (or more) fatty acids. Typically, these fatty acids are randomly distributed among the triglyceride molecules, and the chemist cannot identify anything more than an average composition for a given fat or oil. The average fatty acid composition of some selected fats and oils is given in

Table 1 Common Fatty Acids


C ₁₂ Acids	Lauric	CH ₃ (CH ₂) ₁₀ COOH
C ₁₄ Acids	Myristic	CH ₃ (CH ₂) ₁₂ COOH
C ₁₆ Acids	Palmitic	CH ₃ (CH ₂) ₁₄ COOH
	Palmitoleic	CH ₃ (CH ₂) ₅ CH=CH—CH ₂ (CH ₂) ₆ COOH
C ₁₈ Acids	Stearic	CH ₃ (CH ₂) ₁₆ COOH
	Oleic	CH ₃ (CH ₂) ₇ CH=CH—CH ₂ (CH ₂) ₆ COOH
	Linoleic	CH ₃ (CH ₂) ₄ (CH=CH—CH ₂) ₂ (CH ₂) ₆ COOH
	Linolenic	CH ₃ CH ₂ (CH=CH—CH ₂) ₃ (CH ₂) ₆ COOH
	Ricinoleic	CH ₃ (CH ₂) ₅ CH(OH)CH ₂ CH=CH(CH ₂) ₇ COOH

Table 2. As indicated, all the values in the table may vary in percentage, depending, for instance, on the locale in which the plant was grown or on the particular diet on which the animal subsisted. Thus, perhaps there is a basis for the claims that corn-fed hogs or cattle taste better than animals maintained on other diets.

Vegetable fats and oils are usually found in fruits and seeds and are recovered by three principal methods. In the first method, **cold pressing**, the appropriate part of the dried plant is pressed in a hydraulic press to squeeze out the oil. The second method is **hot pressing**, which is the same as the first method but done at a higher temperature. Of the two methods, cold pressing usually gives a better grade of product (more bland); the hot pressing method gives a higher yield, but with more undesirable constituents (stronger odor and flavor). The third method is **solvent extraction**. Solvent extraction gives the highest recovery of all and can now be regulated to give bland, high-grade food oils.

Animal fats are usually recovered by **rendering**, which involves cooking the fat out of the tissue by heating it to a high temperature. An alternative method involves placing the fatty tissue in boiling water. The fat floats to the surface and is easily recovered. The most common animal fats, lard (from hogs) and tallow (from cattle), can be prepared in either way.

Many triglyceride fats and oils are used for cooking. We use them to fry meats and other foods and to make sandwich spreads. Almost all commercial cooking fats and oils, except lard, are prepared from vegetable sources. Vegetable oils are liquids at room temperature. If the double bonds in a vegetable oil are hydrogenated, the resultant product becomes solid. In making commercial cooking fats (Crisco, Spry, Fluffo, etc.), manufacturers hydrogenate a liquid vegetable oil until the desired degree of consistency is achieved. This makes a product that still has a high degree of unsaturation (double bonds) left. The same technique is used for margarine. "Polyunsaturated" oleomargarine is produced by the partial hydrogenation of oils from corn, cottonseed, peanut, and soybean sources. The final product has a yellow dye (β -carotene) added to make it look like butter; milk, about 15% by volume, is mixed into it to form the final emulsion. Vitamins A and D are also commonly added. Because the final product is tasteless (try Crisco), salt, acetoin, and biacetyl are often added. The latter two additives mimic the characteristic flavor of butter.

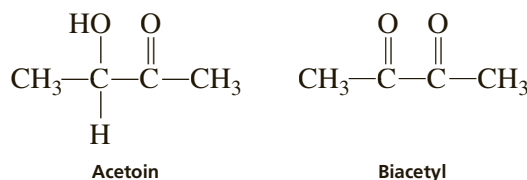

Many producers of margarine claim it to be more beneficial to health because it is "high in polyunsaturates." Animal fats are low in unsaturated fatty acid content and are generally excluded from the diets of people who have high cholesterol levels.


Table 2 Average Fatty Acid Composition (by Percentage) of Selected Fats and Oils

	C ₁₀ C ₈ C ₆ C ₄	C ₁₂ Lauric	C ₁₄ Myristic	C ₁₆ Palmitic	C ₁₈ Stearic	C ₂₀ C ₂₂ C ₂₄	C ₁₆ Palmitoleic	C ₁₈ Oleic	C ₁₈ Ricinoic	C ₁₈ Linoleic	C ₁₈ Linolenic	C ₁₈ Eleostearic	C ₂₀ C ₂₂ C ₂₄		
	Saturated Fatty Acids (No Double Bonds)						Unsaturated (1 Double Bond)						Unsaturated (>1 Double Bond)		
							24-30								
Animal fats															
Tallow			2-3	24-32	14-32		1-3	35-48		2-4					
Butter	7-10	2-3	7-9	23-26	10-13		5	30-40		4-5			2		
Lard			1-2	28-30	12-18		1-3	41-48		6-7			2		
Animal oils															
Neat's foot				17-18	2-3			74-77							
Whale			4-5	11-18	2-4		13-18	33-38					17-31		
Sardine			6-8	10-16	1-2		6-15						12-19		
Vegetable oils															
Corn			0-2	7-11	3-4		0-2	43-49		34-42					
Olive			0-1	5-15	1-4		0-1	69-84		4-12					
Peanut				6-9	2-6	3-10	0-1	50-70		13-26					
Soybean			0-1	6-10	2-6			21-29		50-59	4-8				
Safflower				6-10	1-4			8-18		70-80	2-4				
Castor bean				0-1				0-9	80-92	3-7					
Cottonseed			0-2	19-24	1-2		0-2	23-33		40-48					
Linseed				4-7	2-5			9-38		3-43	25-58				
Coconut	10-22	45-51	17-20	4-10	1-5			2-10		0-2					
Palm			1-3	34-43	3-6			38-40		5-11					
Tung								4-16		0-1			74-91		

Such people have difficulty in metabolizing saturated fats correctly and should avoid them because they encourage cholesterol deposits to form in the arteries. This ultimately leads to high blood pressure and heart trouble. People who pay close attention to their intake of fats tend to avoid consuming large quantities of saturated fats, knowing that eating these fats increases the risk of heart disease. Diet-conscious people try to limit their fat consumption to unsaturated fats, and they make use of the current mandatory food labeling to obtain information on the fat content of the food they eat.

Unfortunately, not all of the unsaturated fats appear to be equally safe. When we eat partially hydrogenated fats, we increase our consumption of *trans*-fatty acids. These acids, which are isomers of the naturally occurring *cis*-fatty acids, have been implicated in a variety of conditions, including heart disease, cancer, and diabetes. The strongest evidence that *trans*-fatty acids may be harmful comes in studies of the incidence of coronary heart disease. Ingestion of *trans*-fatty acids appears to increase blood cholesterol levels, in particular the ratio of low-density lipoproteins (LDL, or “bad” cholesterol) to high-density lipoproteins (HDL, or “good” cholesterol). The *trans*-fatty acids appear to exhibit harmful effects on the heart that are similar to those shown by saturated fatty acids.


The *trans*-fatty acids do not occur naturally to any significant extent. Rather, they are formed during the partial hydrogenation of vegetable oils to make margarine and solid forms of shortening. For a small percentage of *cis*-fatty acids subjected to hydrogenation, only one hydrogen atom is added to the carbon chain. This process forms an intermediate free radical, which is able to rotate its conformation by 180 degrees before it releases the extra hydrogen atom back to the reaction medium. The result is an isomerization of the double bond.

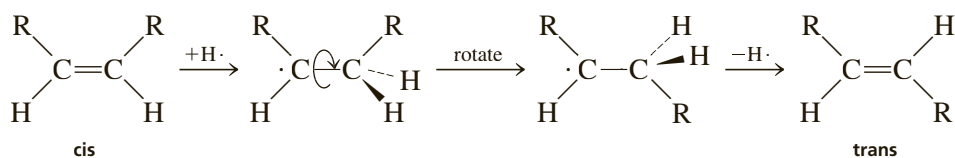


Concern over the health and nutrition of the public, particularly over the average fat intake of most Americans, has prompted food chemists and technologists to develop a variety of **fat replacers**. The objective has been to discover substances that have the taste and mouth-feel of a real fat, but do not have deleterious effects on the cardiovascular system. One product that has been used in certain snack foods is **olestra** (marketed under the trade name **Olean**, by the Procter and Gamble Company). Olestra is not an acylglycerol; rather, it is composed of a molecule of **sucrose** that has been substituted by long-chain fatty acid residues. It is a **polyester**, and the body's enzyme systems are not capable of attacking it and catalyzing its breakdown into smaller molecules.

Because the body's enzyme systems are unable to break this molecule down, it does not contain any usable dietary calories. Furthermore, it is heat stable, which makes it ideal for frying and other cooking. Unfortunately, for some individuals there may be harmful or unpleasant side effects. The use of olestra has been reported to deplete certain fat-soluble vitamins, particularly vitamins A, D, E, and K. For this reason, products prepared with olestra have these vitamins added to offset this effect. Also, some people have reported diarrhea and abdominal cramps. Currently, only a few products that contain olestra are available in the U.S.

Is the development of fat replacers such as olestra part of the wave of the future? As the average American's appetite for snack foods continues to grow and as health problems arising from obesity also increase, the demand for satisfying foods that are less fattening will always be strong. In the long run, however, it would probably be better if we all learned to curtail our appetite for fatty foods and, instead, tried to increase our intake of fruits, vegetables, and other healthful foods.

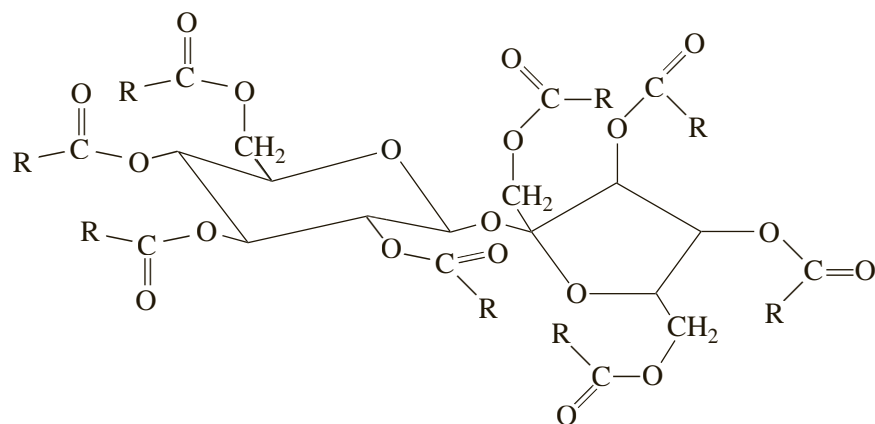

Olestra

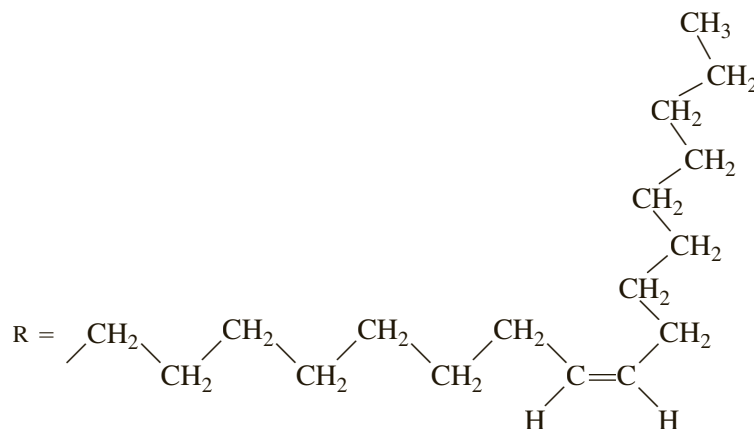

At the same time, a change from a sedentary lifestyle to one that include regular exercise would also be much more beneficial to our health.

REFERENCES

- Dawkins, M. J. R.; Hull, D. The Production of Heat by Fat. *Sci. Am.* **1965**, 213 (Aug), 62.
- Dolye, E. Olestra? The Jury's Still Out. *J. Chem. Educ.* **1997**, 74 (Apr), 370.
- Eckey, E. W.; Miller, L. P. *Vegetable Fats and Oils*; ACS Monograph 123; Reinhold: New York, 1954.
- Farines, M.; Soulier, F.; Soulier, J. Analysis of the Triglycerides of Some Vegetable Oils. *J. Chem. Educ.* **1988**, 65 (May), 464.
- Gunstone, F. D. The Composition of Hydrogenated Fats Determined by High Resolution ¹³C NMR Spectroscopy. *Chem. Ind.* **1991**, (Nov 4), 802.
- Heinzen, H.; Moyna, P.; Grompone, A. Gas Chromatographic Determination of Fatty Acid Compositions. *J. Chem. Educ.* **1985**, 62 (May), 449.
- Jandacek, R. J. The Development of Olestra, a Noncaloric Substitute for Dietary Fat. *J. Chem. Educ.* **1991**, 68 (Jun), 476.
- Kalbus, G. E.; Lieu, V. T. Dietary Fat and Health: An Experiment on the Determination of Iodine Number of Fats and Oils by Coulometric Titration. *J. Chem. Educ.* **1991**, 68 (Jan), 64.
- Lemonick, M. D. Are We Ready for Fat-Free Fat? *Time* **1996**, 147 (Jan 8), 52.
- Martin, C. TFA's—a Fat Lot of Good? *Chem. Brit.* **1996**, 32 (Oct), 34.
- Nawar, W. W. Chemical Changes in Lipids Produced by Thermal Processing. *J. Chem. Educ.* **1984**, 61 (Apr), 299.
- Shreve, R. N.; Brink, J. Oils, Fats, and Waxes. *The Chemical Process Industries*, 4th ed.; McGraw-Hill: New York, 1977.
- Thayer, A. M. Food Additives. *Chem. Eng. News* **1992**, 70 (Jun 15), 26.
- Wootan, M.; Liebman, B.; Rosofsky, W. *trans*: The Phantom Fat. *Nutr. Action Health Lett.* **1996**, 23 (Sep), 10.


Methyl Stearate from Methyl Oleate

Catalytic hydrogenation

Recrystallization


Unsaturation tests

In this experiment, you will convert the liquid methyl oleate, an “unsaturated” fatty acid ester, to solid methyl stearate, a “saturated” fatty acid ester, by catalytic hydrogenation.

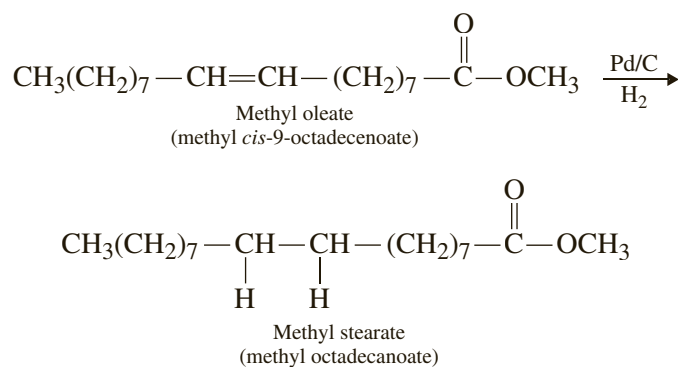

By commercial methods such as those described in this experiment, the unsaturated fatty acids of vegetable oils are converted to margarine (see the essay “Fats and Oils”). However, rather than using the mixture of triglycerides that would be present in a cooking oil such as Mazola (corn oil), we use as a model the pure chemical methyl oleate.

For this procedure, a chemist would usually use a cylinder of hydrogen gas. Because many students will be following the procedure simultaneously, however, we use the simpler expedient of causing zinc metal to react with dilute sulfuric acid:

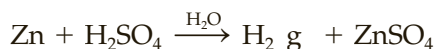

The hydrogen so generated will be passed into a solution containing methyl oleate and the palladium on carbon catalyst (10% Pd/C).

REQUIRED READING

Review: Techniques 1–6

New: Technique 8 Filtration, Sections 8.3–8.5

Technique 9 Physical Constants of Solids: The Melting Point

Essay Fats and Oils

You should also read those sections in your lecture textbook that deal with catalytic hydrogenation. If the instructor indicates that you should perform the optional

unsaturation tests on your starting material and product, read the descriptions of the $\text{Br}_2/\text{CH}_2\text{Cl}_2$ test in Experiment 57C and in Experiment 24.

SPECIAL INSTRUCTIONS

Because this experiment calls for generating hydrogen gas, no flames will be allowed in the laboratory.

CAUTION



No flames allowed.

Because a buildup of hydrogen is possible within the apparatus, it is especially important to remember to wear your safety goggles; you can thus protect yourself against the possibility of minor “explosions” from joints popping open, from fires, or from any glassware accidentally cracking under pressure.

CAUTION



Wear safety goggles.

When you operate the hydrogen generator, be sure to add sulfuric acid at a rate that does not cause hydrogen gas to evolve too rapidly. The hydrogen pressure in the vial should not rise much above atmospheric pressure. Neither should the hydrogen evolution be allowed to stop. If this happens, your reaction mixture may be “sucked back” into your hydrogen generator.

SUGGESTED WASTE DISPOSAL

Dispose of the sulfuric acid (from the hydrogen generator) by pouring it into the waste container designated for acid waste. Place any leftover zinc in the container designated for this material. After centrifugation, transfer the Pd/C catalyst to a container designated for this material. After collecting the methyl stearate by filtration, place the methanol filtrate in the nonhalogenated organic waste container. Discard the solutions that remain after the bromine test for unsaturation into a waste container designated for the disposal of halogenated organic solvents.

NOTES TO THE INSTRUCTOR

Use methyl oleate that is 100% or nearly 100% pure. We use Aldrich Chemical Co. No. 311111.

PROCEDURE

Apparatus


Assemble the apparatus as illustrated in the figure. The apparatus can be simplified by using the multipurpose adapter (Figure 14.9) in place of the Claisen

head and both thermometer adapters. The apparatus consists of basically three parts:

1. Hydrogen generator
2. Reaction flask
3. Mineral oil bubbler trap

The mineral oil bubbler trap has two functions. First, it allows you to keep a pressure of hydrogen within the system that is slightly above atmospheric. Second, it prevents back-diffusion of air into the system. The functions of the other two units are self-explanatory.

So that hydrogen leakage is prevented, the tubing used to connect the various subunits of the apparatus should be either relatively new rubber tubing, without cracks or breaks, or Tygon tubing. The tubing can be checked for cracks or breaks simply by stretching and bending it before use. It should be of such a size that it will fit onto all connections tightly. Similarly, if any rubber stoppers are used, they should be fitted with a size of glass tubing that fits firmly through the holes in their centers. If the seal is tight, it will not be easy to slide the glass tubing up and down in the hole. The inlet tube (Pasteur pipette) in the round-bottom flask should reach almost to the bottom of the flask. Hydrogen must bubble *through* the solution.

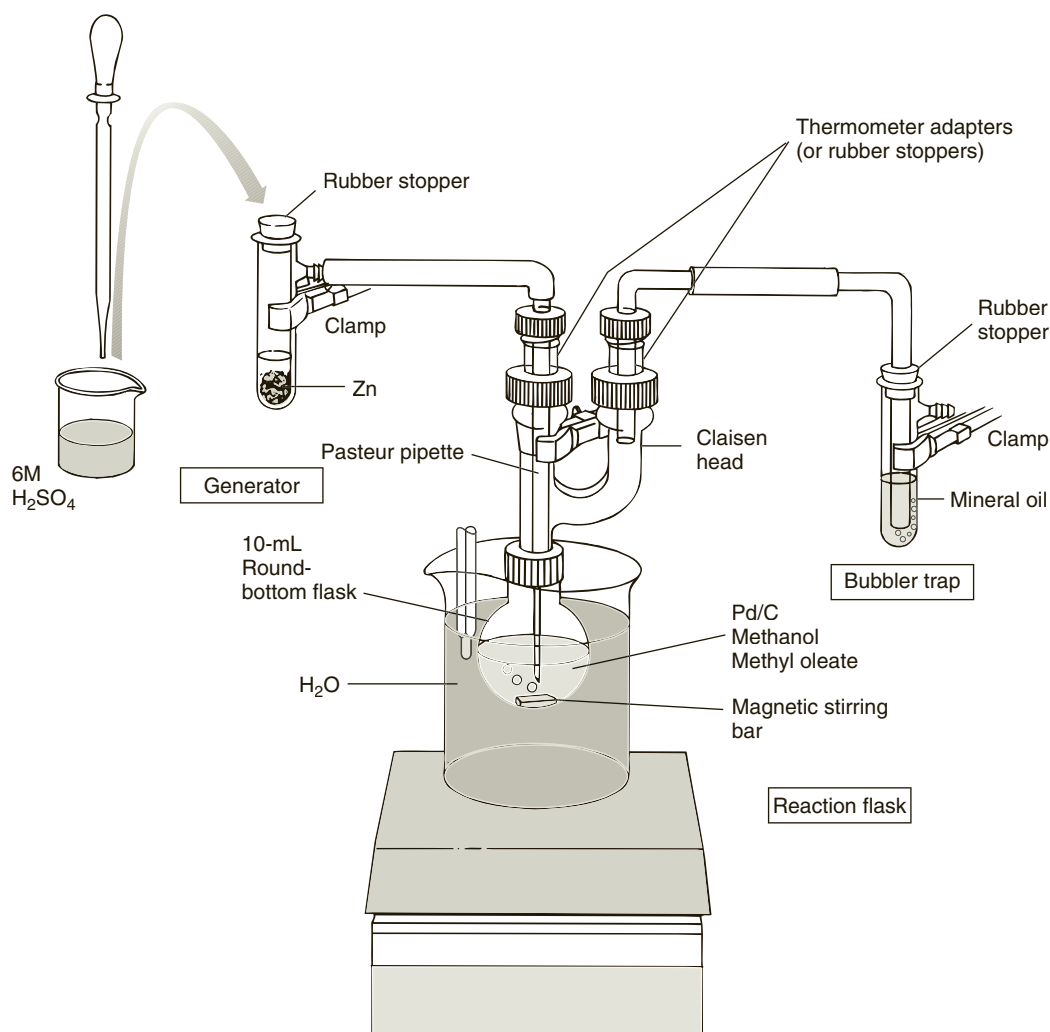

Apparatus for Experiment 25.

Preparing for the Reaction

Fill the bubbler trap (second side-arm test tube) about one-third full with mineral oil. The end of the glass tube should be submerged below the surface of the oil.

To charge the hydrogen generator, weigh out about 2 g of mossy zinc and place it in the side-arm test tube. Seal the large opening at its top using a rubber stopper. Obtain about 10 mL of 6M sulfuric acid and place it in a small Erlenmeyer flask or beaker **but do not add it yet**.

Weigh a 10-mL round-bottom flask and then place 1.00 mL of methyl oleate into it. Reweigh the flask to obtain the exact amount of methyl oleate used. After this, add 6.0 mL of methanol solvent to the flask. Also place a magnetic stirring bar into the flask. Place the flask into a small beaker. Using smooth weighing paper, weigh about 0.030 g (30 mg) of 10% Pd/C. Carefully place about one-third of the catalyst into the flask and gently swirl the liquid until the solid catalyst has sunk into the liquid. Repeat this with the rest of the catalyst, adding one-third of the original amount each time.

CAUTION



Be careful when adding the catalyst; sometimes it will cause a flame. Do not hold on to the flask; it should be in a small beaker on the lab bench. Have a watch glass handy to cover the opening and smother the flame should this occur.

Running the Reaction

Complete the assembly of the apparatus, making sure that all the seals are gas tight. Place the round-bottom flask in a warm-water bath maintained at 40°C. This will help to keep the product dissolved in the solution throughout the course of the reaction. If the temperature rises above 40°C, you will lose a significant amount of the methanol solvent. If this occurs, do not hesitate to add more methanol to the reaction flask through the side arm of the Claisen head, using a Pasteur pipette. Begin stirring the reaction mixture with the magnetic stirring bar. Start the evolution of hydrogen by removing the rubber stopper and adding a portion of the 6M sulfuric acid solution (about 6 mL) to the hydrogen generator. Replace the rubber stopper. A good rate of bubbling in the reaction flask is about three to four bubbles a second. Continue the evolution of hydrogen for 45–60 minutes. If necessary, open the generator, empty it, and refresh the zinc and sulfuric acid. (Keep in mind that the acid is used up and becomes more dilute as the zinc reacts.)

Stopping the Reaction

After the reaction is complete, stop the reaction by disconnecting the reaction flask. Decant the acid in the side-arm test tube into a designated waste container, being careful not to transfer any zinc metal. Rinse the zinc in the test tube several times with water and then place any unreacted zinc in a waste container provided for this purpose.

Keep the temperature of the reaction mixture at about 40°C until you perform the centrifugation; otherwise, the methyl stearate may crystallize and interfere with removal of the catalyst. There should not be any white solid (product) in the round-bottom flask. If there is a white solid, add more methanol and stir until the solid dissolves.

Removal of the Catalyst

Pour the reaction mixture into a centrifuge tube. Place the centrifuge tube into the water bath at 40°C until just before you are ready to centrifuge the mixture. Centrifuge the mixture for several minutes. After centrifugation, the black catalyst should

be at the bottom of the tube. If some of the catalyst is still suspended in the liquid, heat the mixture to 40°C and centrifuge the mixture again. Carefully pour the supernatant liquid (leaving the black catalyst in the centrifuge tube) into a small beaker and cool to room temperature.

Crystallization and Isolation of Product

Place the beaker in an ice bath to induce crystallization. If crystals do not form or if only a few crystals form, you may need to reduce the volume of solvent. Do this by heating the beaker in a water bath and directing a slow stream of air into the beaker, using a Pasteur pipette for a nozzle (Figure 7.12A). If crystals begin to form while you are evaporating the solvent, remove the beaker from the water bath. If crystals do not form, reduce the volume of the solvent by about one third. Allow the solution to cool and then place it in an ice bath.

Collect the crystals by vacuum filtration, using a small Hirsch funnel (Technique 8, Section 8.3). Save both the crystals and the filtrate for the tests later. After the crystals are dry, weigh them and determine their melting point (literature value, 39°C). Calculate the percentage yield. Submit your remaining sample to your instructor in a properly labeled container, along with your report.

Unsaturation Tests (Optional)

Number three test tubes (1, 2, 3) and place one of the following samples into each test tube:

1. About 0.1 mL of methyl oleate dissolved in a small amount of methylene chloride
2. A small spatulaful of your methyl stearate product dissolved in a small amount of methylene chloride
3. About 0.1 mL of the filtrate that you saved as directed earlier


To each test tube, add a solution of bromine in methylene chloride, drop by drop, to the contents of the test tubes until the red color is no longer discharged. Record the results in each case, including the number of drops required. What do these results indicate about the presence of unsaturated compounds in each sample?

QUESTIONS

1. Using the information in the essay on fats and oils, draw the structure of the triacylglycerol (triglyceride) formed from oleic acid, linoleic acid, and stearic acid. Give a balanced equation and show how much hydrogen would be needed to reduce the triacylglycerol completely; show the product.
2. A 0.150-g sample of a pure compound subjected to catalytic hydrogenation takes up 25.0 mL of H₂ at 25°C and 1 atm pressure. Calculate the molecular weight of the compound, assuming that it has only one double bond.
3. A compound with the formula C₅H₆ takes up 2 moles of H₂ on catalytic hydrogenation. Give one possible structure that would fit the information given.
4. A compound of formula C₆H₁₀ takes up 1 mole of H₂ on reduction. Give one possible structure that would fit the information.


Soap

Soaps as we know them today were virtually unknown before the first century AD. Clothes were cleaned primarily by the abrasive action of rubbing them on rocks in water. Later, it was discovered that certain types of leaves, root, nuts, berries, and barks formed soapy lathers that solubilized and removed dirt from clothes. We now refer to these natural materials that lather as **saponins**. Many saponins contain pentacyclic triterpene carboxylic acids, such as oleanolic acid or ursolic acid, chemically combined with a sugar molecule. These acids also appear in the uncombined state. Saponins were probably the first known “soaps.” They may have also been an early source of pollution in that they are known to be toxic to fish. The pollution problem associated with the development of soap and detergents has been long and controversial.

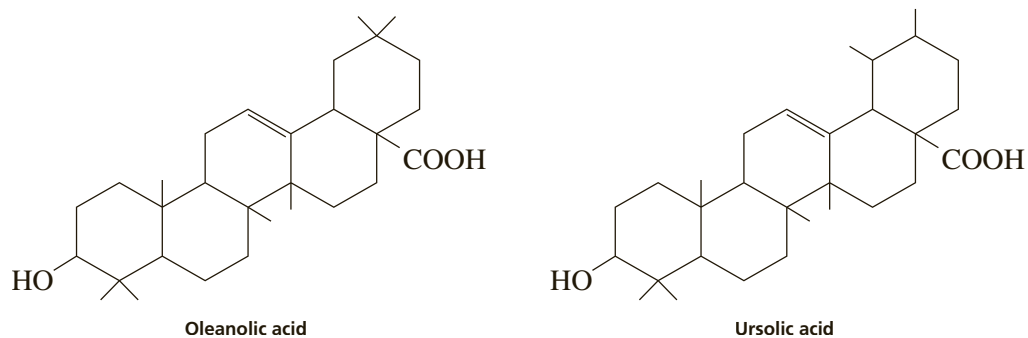

Soaps today have evolved over many centuries from experimentation with crude mixtures of alkaline and fatty materials. Pliny the Elder described the manufacture of soap during the first century AD. A modest soap factory was even built in Pompeii. During the Middle Ages, cleanliness of the body or clothing was not considered important. Those who could afford perfumes used them to hide their body odor. Perfumes, like fancy clothes, were status symbols for the rich. An interest in cleanliness again emerged during the eighteenth century, when disease-causing microorganisms were discovered.


The process of making soap has remained practically unchanged for 2000 years. The procedure involves the basic hydrolysis or **saponification** of an animal fat or a vegetable oil. Chemically, fats and oils are referred to as **triglycerides** or **triacylglycerols**. They contain ester functional groups. Saponification involves heating a fat or oil with an alkaline solution. This alkaline solution was originally obtained by the leaching of wood ashes or from the evaporation of natural alkaline waters. Today, lye (sodium hydroxide) is used as the source of the alkali. The alkaline solution hydrolyzes the fat or oil into its component parts, the sodium salt of a long-chain carboxylic acid (soap) and an alcohol (glycerol). The following equation shows how soap is produced from a fat or oil.

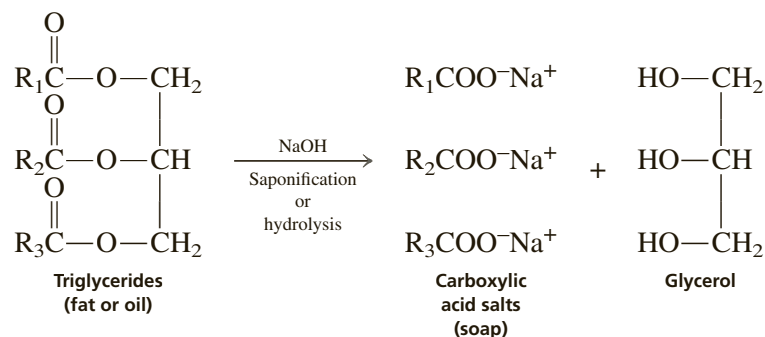

The fatty acids are rarely of a single type in any given fat or oil. In fact, a single triglyceride molecule in a fat may contain three different acid residues (R_1COOH , R_2COOH , R_3COOH), and not every triglyceride in the substance will be identical. Each fat or oil, however, has a characteristic *statistical distribution* of the various types of acids possible. The carboxylic acid salts of soap usually contain 12–18 carbons arranged in a straight chain. Soaps made with short-chain fatty acids are irritating to skin, while saturated fatty acids with long chain lengths are not as soluble. The carboxylic acids containing even numbers of carbon atoms predominate, and the chains may contain unsaturation. The composition of the common fats and oils is given in the essay “Fats and Oils,” which precedes Experiment 23.

The fats and oils that are most common in commercial soap preparations are lard and tallow from animal sources, and coconut, palm, and soybean oil from vegetable sources. The length of the hydrocarbon chain and the number of double bonds in the carboxylic acid portion of the fat or oil determine the properties of the resulting soap.

Soaps can be made from a single fat or oil, but it is much more common to use two or more fats or oils in order to make soap with desired characteristics. For example, olive oil is widely used by home and small-scale soap makers. 100% olive oil soap (castile soap) is very moisturizing due to its high oleic acid content, but the bar is softer, has a shorter shelf and use life and can leave behind an oily feeling. Coconut oil, having a high lauric acid content, produces a hard, bubbly bar, with a longer shelf life with excellent cleansing and antimicrobial properties. However, bars with more than 40% coconut oil can be drying to the skin. Therefore, a soap made from 30% coconut oil and 70% olive oil will exhibit cleaning quality of coconut oil as well as the moisturizing quality of olive oil, producing a more desirable soap than soap prepared from 100% of either one of these oils. All of the soap oil combinations in Experiment 24 contain coconut oil as it is key to making a good soap.

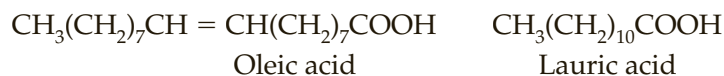

Tallow and coconut oil are the principal fatty materials used in making soap. The cattle industry makes tallow cheap and readily available. It is well suited for soap due to its long-chain saturated and unsaturated fatty acid content, adding hardness and creaminess to the soap. Rendering tallow from the solid fat of cattle (suet) can be as simple as chopping, heating, and straining. Commercial processes often involve melting with steam and multiple stages of centrifugal separation.


Commercial and home soap makers blend tallow with coconut oil and saponify this mixture. A similar soap can be made from lard (from hogs), which differs from tallow in that lard contains more oleic acid. In Experiment 24A, you will make a soap consisting of 70% lard and 30% coconut oil.

When soap is prepared from an oil or fat and base, glycerol is also produced. Glycerol is a humectant, which attracts moisture from the air, and an emollient, which softens the skin. Because of these desirable characteristics, commercial saponification will often remove the glycerol for use in the cosmetics industry, while home soap makers and small-scale soap companies will usually leave the glycerol in the soap. To give soap a desirable scent, essential oils or perfumes are often added to the soap before it is poured into a mold.


Saponification to make soap can take place either with sodium hydroxide or potassium hydroxide. The molecules crystallize differently depending on the base used. Sodium hydroxide forms a tighter crystal matrix, producing a harder bar while potassium hydroxide forms more loosely, and is used more frequently for liquid soaps.

The two most commonly used methods to make soap at home are called the **cold process** and the **hot process**. Both require a heat source and careful calculations to ensure that no caustic base is left unreacted in the soap. The hot process uses more heat to speed the reaction, resulting in a fully saponified soap by the time the soap is put into molds. The cold process uses just enough heat to ensure that all the fat is melted prior to reacting it with the base. We will use the cold process in Experiment 24 because it is simpler, requires less laboratory time, and results in a creamier bar. The disadvantage of the cold process is that complete saponification requires from two days to a week or longer. If the soap is used before complete saponification has occurred, it may cause skin irritation or damage because of the high pH of unreacted base.

A disadvantage of soap is that it is an ineffective cleanser in hard water. Hard water contains salts of magnesium, calcium, and iron in solution. When soap is used in hard water, “calcium soap,” the insoluble calcium salts of the fatty acids, and other precipitates are deposited as curds. This precipitate, or curd, is referred to as bathtub ring. Although soap is a poor cleanser in hard water, it is an excellent cleanser in soft water.

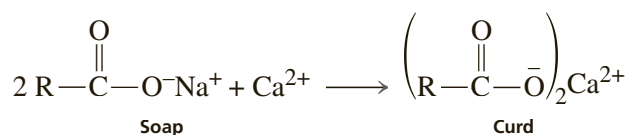


Water softeners are sometimes added to soaps to help remove the troublesome hard-water ions so that the soap will remain effective in hard water. Sodium carbonate or trisodium phosphate will precipitate the ions as the carbonate or phosphate salts. Unfortunately, the precipitate may become lodged in the fabric of items being laundered, causing a grayish or streaked appearance.

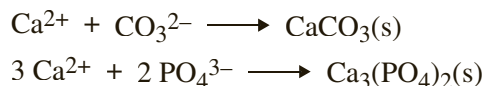

An important advantage of soap is that it is **biodegradable**. Microorganisms can consume the linear soap molecules, converting them to carbon dioxide and water, eliminating soap from the environment.

ACTION OF SOAP IN CLEANING

Dirty clothes, skin, or other surfaces have particles of dirt suspended in a layer of oil or grease. Polar water molecules cannot remove the dirt embedded in nonpolar oil or grease. One can remove the dirt with soap, however, because of its dual nature. The soap molecule is a surfactant, having a polar, *water-soluble* head (carboxylate salt) and a long, *oil-soluble* tail (hydrocarbon chain). The hydrocarbon tail of soap dissolves in the oily substance, but the ionic end remains outside the oily surface. Soap molecules orient themselves such that their hydrocarbon ends dissolve into the soiled oil. The negatively charged oil droplet becomes strongly attracted to, and solvated by, water, which lifts and carries the soiled oil away from the surface. The solvated oil droplet is called a **micelle**.

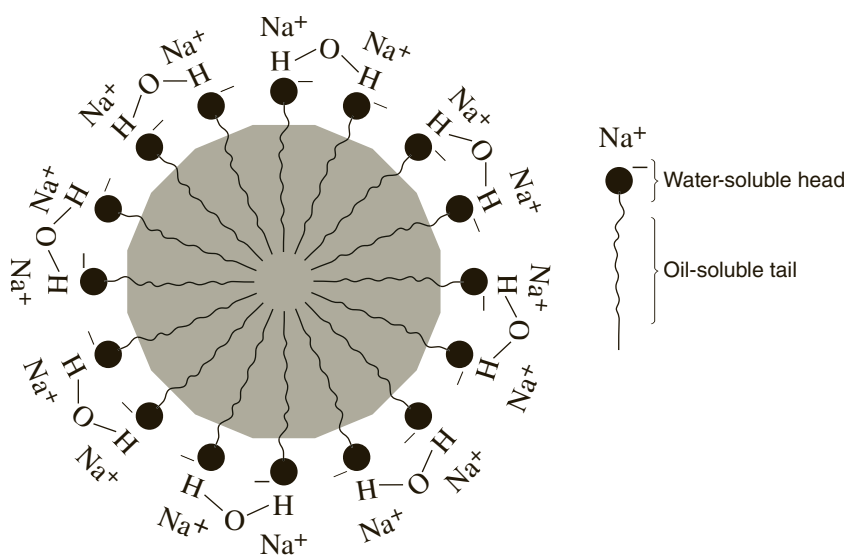


A Soap micelle solvating a droplet of oil (from Linstromberg, W. W., *Organic Chemistry: A Brief Course*, D. C. Heath, 1978).

REFERENCES

- Ainsworth, S. J. Soaps and Detergents. *Chem. Eng. News* **1994**, 72 (Jan 24), 34.
 Ainsworth, S. J. Soaps and Detergents. *Chem. Eng. News* **1995**, 73 (Jan 23), 30.
 Ainsworth, S. J. Soaps and Detergents. *Chem. Eng. News* **1996**, 74 (Jan 23), 34.
 Snell, F. D. Soap and Glycerol. *J. Chem. Educ.* **1942**, 19 (Apr), 172.
 Thayer, A. M. Soaps and Detergents. *Chem. Eng. News* **1993**, 71 (Jan 25), 26.

Preparation of Soap¹

In this experiment we prepare soap from animal fats and vegetable oils. Fats and oils (sometimes called **triglycerides**) are esters derived from three fatty (carboxylic) acids and glycerol. Base hydrolysis of triglycerides (**saponification**) cleaves the ester bond, releasing the fatty acid salts (soap) and glycerol. (See the preceding essay for a complete discussion of soaps.)

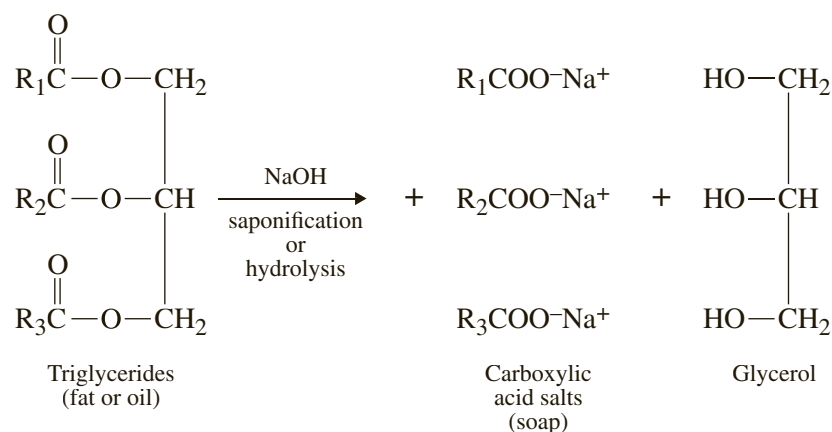

To calculate the amount of sodium hydroxide required to make soap with a specific composition, you must know the **saponification value** of each fat or oil in the soap. The saponification value or number (SAP) is defined as the number of milligrams of potassium hydroxide required to saponify one gram of fat or oil. Saponification values can be converted to number of grams of sodium hydroxide required to saponify one gram of fat or oil (see table below).

Fat or oil	Amount of NaOH (g) required to saponify one gram of the fat or oil
Avocado oil	0.133
Canola oil	0.132
Coconut oil	0.181
Lard	0.141
Olive oil	0.135
Vegetable shortening (Crisco)	0.136

For example, 0.133 g of solid sodium hydroxide would fully saponify 1.0 g of avocado oil. In this experiment, you will use a 5% excess of fat or oil to ensure all of the sodium hydroxide is consumed in the saponification process. Unreacted

¹This experiment was developed by Sean Ichun Choi, organic chemistry student at North Seattle Community College and soap maker extraordinaire!

sodium hydroxide increases the pH and could produce soap harmful to human skin. The excess unreacted fat or oil is referred to as the **superfat percentage** or **lye discount**. If 100 g of avocado oil were used, then the amount of sodium hydroxide required (assuming a 5% superfat) would be:

$$95 \text{ g oil} \times (0.133 \text{ g NaOH}/1.0 \text{ g oil}) = 12.63 \text{ g of solid NaOH.}$$

Water is also required to dissolve the sodium hydroxide and to ensure thorough mixing of sodium hydroxide and fat or oil. In this experiment, the amount of water to be used will be determined by multiplying the total weight of sodium hydroxide by 2. For a soap made with 100 g of avocado oil, the amount of water required would be $2 \times 12.63 \text{ g} = 25.3 \text{ g}$ (25.3 mL of water).

Calculating the amount of sodium hydroxide to use when making soap comprised of two or more different fats or oils is slightly more complicated. Each fat or oil needs to be calculated individually using their respective SAP values (see table) and then added together. Here are the calculations for the soap that will be made in Experiment 26A, 30% coconut oil and 70% lard. For these calculations, we assume that 30 g of coconut oil and 70 g of lard will be used. First, calculate the mass of 95% of each oil or fat to ensure that there is a 5% excess of oil (5% superfat). Then, do the following calculation:

$$[28.5 \text{ g coconut oil} \times 0.181 \text{ g NaOH/g oil}] + [66.5 \text{ g lard} \times 0.141 \text{ g NaOH/g oil}] \\ = 5.16 \text{ g} + 9.38 \text{ g} = 14.54 \text{ g NaOH required.}$$

The amount of water to be used is determined by multiplying the total weight of sodium hydroxide by two:

$$2 \times 14.54 \text{ g} = 29.1 \text{ g (29.1 mL) water required}$$

The recipe would be 30 g of coconut oil, 70 g of lard, 14.54 g sodium hydroxide, and 29.1 mL of water. Remember that although the calculations are performed using 95% of the fat or oil weight, the recipe calls for 100% of the weight of the fat or oil. In Experiment 26A, the amount of soap prepared will be reduced by dividing all of these values by 3.

In Experiment 26B, the compositions of several different soaps are presented as options. You will need to calculate how much sodium hydroxide and water are required for a given soap using the saponification table.

REQUIRED READING

New: Essay Soap

SPECIAL INSTRUCTIONS

It is very important that you check the pH of your soap before removing it from the mold and before taking the soap home or using it. The best way to measure the pH is with a meter. Soaps with $\text{pH} > 10$ should not be used.

NOTES TO THE INSTRUCTOR

It is best to use ground sodium hydroxide rather than pellets. This makes it much easier for students to weigh out the calculated amount of sodium hydroxide. This should be done quickly to avoid absorption of water.

Students should wait at least 2 or 3 weeks before removing their soap from the mold. It is very important that students check the pH of their soap before removing it from the mold and before taking their soap home or using it. The best way to measure the pH is with a meter. Many pH papers such as Hydrion are unreliable with soap solutions because the surfactant nature of soap interacts with the water/dye in the pH paper. This usually results in a pH reading lower than the actual value. ColorpHast® pH Test Strips (#9580) are more expensive, but provide more accurate pH readings with soap, and are suitable for this experiment. Phenolphthalein is an economical way to test pH more reliably than pH paper. Because the phenolphthalein test does not give a precise pH, students should also test the pH either with a meter or the ColorpHast® pH Test Strips. Soaps with pH > 10 should not be used.

Weighing boats provide a convenient mold for the soap. We use weighing boats from Cynmar, Stock # 11120957, 3.5 in × 3.5 in × 1 in. To make it easier to remove the bar of soap, apply a thin coating of petroleum jelly to the inside of the weighing boat before pouring the soap mixture into the weighing boats. If students use other containers for molds, the containers should hold at least 50 mL of liquid. Do **not** use an aluminum mold or liner since it will react with the sodium hydroxide to form flammable hydrogen gas.

SUGGESTED WASTE DISPOSAL

Any soap with a pH greater than 10 should be discarded into a container specifically designated for this material. Dispose of any left over soap solutions (from the pH or other tests) in the container for aqueous waste.

26A EXPERIMENT 26A

Preparation of Soap from 70% Lard and 30% Coconut Oil

In this experiment you will prepare soap using a procedure that gives the specific quantities of all required materials. The soap prepared in this experiment is a basic, all-purpose soap.

PROCEDURE

CAUTION



Sodium hydroxide may be supplied as pellets or granules. Do not touch the pellets or granules since sodium hydroxide is very caustic and can cause chemical burns. Use a scoop or forceps and wear gloves. Sodium hydroxide should be weighed as rapidly as possible since it is hygroscopic and draws moisture from the air, becoming sticky. Keep the lid on the bottle as much of the time as possible. If the sodium hydroxide is in pellet form, you will not be able to weigh out exactly 4.85 g. Don't try to break up a pellet—just get as close to 4.85 g as possible with whole pellets. When making up the sodium hydroxide solution, you should always add sodium hydroxide to water. Adding water to sodium hydroxide can cause an intense localized temperature increase that can boil the liquid and cause splattering.

Preparation of Reagents. Add 9.7 mL of distilled water to a 50-mL Erlenmeyer flask along with a magnetic stir bar. While stirring, slowly add 4.85 g of solid sodium hydroxide to the flask. Continue to stir until all of the sodium hydroxide is dissolved. Note that considerable heat is released during this step. Before adding to the lard and coconut oil, the sodium hydroxide solution should be allowed to cool to 45°C–50°C. If the solution cools below 45°C, you should heat the solution on a hot plate to achieve a temperature within 45°C–50°C.

Add 23.3 g of lard and 10.0 g of coconut oil to a 150-mL beaker along with a magnetic stirring bar. Heat and stir this mixture until all of the solids have melted and the temperature is 45°C–50°C.

Saponification. When both the sodium hydroxide solution and the melted fat and oil mixture are 45°C–50°C, turn off any heat source and slowly add the aqueous base solution with a Pasteur pipette to the fat–oil mixture while stirring rapidly with the magnetic stirrer. Rapid stirring is necessary to thoroughly mix all the ingredients into one layer so they will react more completely. Continue to stir the mixture rapidly until the mixture has thickened to where you can see “trails” of your mixture on the surface. This stage of the saponification process is called **trace**. Trace (trailing) can be observed by dipping a stirring rod or spatula into the mixture, then allowing drops of the mixture to fall back into the beaker as you pull the stirring rod out of the liquid. Look for a small streak that is slightly elevated relative to the surface of the mixture. If you don’t observe trailing, or if it is not obvious that the mixture has thickened, just stir rapidly for 60 minutes and before moving on to the next step.

Obtain a mold for your soap, most likely a weighing boat. Applying a thin coat of petroleum jelly to the inside wall will make it easier to remove the soap. Label the mold on the bottom with your name. If you use a different kind of mold, your instructor will give you instructions on how to prepare it. Any mold used should have a capacity of about 50 mL. Do **not** use an aluminum mold or liner, since it will react with sodium hydroxide to form flammable hydrogen gas.

Once your mixture has reached “trace,” has thickened considerably, or has reacted for 60 minutes, pour the mixture into the mold. Be mindful not to pour your stir bar into the mold. If you do this, use a second magnetic stir bar or forceps to pull out the first one.

NOTE: A helpful hint for cleaning the beaker is to let it sit until next laboratory period. By that time, the beaker will be lined with soap rather than an oily caustic film!

After about 1 hour, your mixture should have firmed up enough so that the mold can be moved easily without spilling the mixture. Any soap at this stage is caustic due to the unreacted sodium hydroxide. Moving molds filled with less viscous soap mixtures increases the chance of spillage. Always wear gloves and goggles, and handle with caution. The mold should be placed in your lab drawer or in a designated location as directed by your laboratory instructor.

CAUTION



At this point, saponification has not been completed, and the mixture will still be strongly basic. Do not try to remove the soap from the mold and do not take your soap home.

Your soap should be almost completely saponified in 24–48 hours. You may see a whitish powder-like substance on the surface of your soap caused by surface aqueous sodium hydroxide reacting with carbon dioxide in the air to form sodium carbonate. Although it may look unsightly, it is harmless.

Before removing the bar of soap from the mold or taking your soap home, you should check the pH as described in the next section. It is also best to let the soap cure for 4–6 weeks to allow trace amounts of remaining base to react, and excess water to evaporate, resulting in a smoother, longer lasting bar of soap.

pH Test. Before using your soap or taking your soap home, you must test the pH to make sure that it is not too basic. Test the pH of the soap after it has set in the mold for about 1 hour and again after 1 week, 2 weeks, and 3 weeks. If the pH is still high (>10) after 3 weeks, discard the soap since it could cause harm to your body. You should include the following table with your entries for the pH values in your report.

Time Interval	Initial pH	1 week	2 weeks	3 weeks
pH Value				

The most accurate way to test the pH is with a pH meter. Using a spatula, remove about 0.5 g of soap in several small pieces from the bar of soap and add to 10 mL of water in a small beaker. Stir thoroughly until the soap is mostly dissolved, then dip the pH electrode into the mixture and record the pH. To be safe, the pH of your soap should be between 8 and 10. **Soaps with a pH in the 10–11 range can be irritating to the skin and should not be used. Soaps with a pH > 11 can cause damage to the body and should not be used.**

Alternatively, the pH can be tested with a ColorpHast[®] pH test strip. Dip the entire strip into the soap solution until the colors are no longer changing (usually a few seconds). To determine the pH, match up the four colors on the strip with the color standards provided with the strips.

The pH can also be checked with phenolphthalein. Using the same mixture described above, add a couple drops of phenolphthalein solution and stir. The test range for phenolphthalein is about pH 8–10, changing from colorless to a deep pink. Colorless to light pink is fine. The more intense the pink, the more basic it is. A soap that gives a deep pink color could be dangerous to use, and you should verify the pH with a pH meter or a ColorpHast[®] pH test strip.

Tests on Soaps. Place about 0.05 g of your soap in a clean 10-mL graduated cylinder. Add 3 mL of deionized water. Cover the opening with parafilm. Place your thumb over the opening and shake the mixture vigorously for about 15 seconds. After about 30 seconds, observe the level of the foam. Add two drops of 4% calcium chloride solution to the soap mixture from a Pasteur pipette. Cover the opening with parafilm and shake the mixture for 15 seconds. Allow it to stand for 30 seconds. Observe the effect of the calcium chloride on the foam. Do you observe anything else? Add 0.5 g of trisodium phosphate and shake the mixture again for 15 seconds. After 30 seconds, again observe the results. Explain the results of these tests in your notebook.

26B EXPERIMENT 26B

Preparation of Several Soaps with a Given % Composition

In this experiment you will prepare a soap with a given fat–oil composition. You will not be given the specific quantities of all required materials. Instead you must calculate the quantities using the table given in Experiment 26A.

PROCEDURE

CAUTION



Sodium hydroxide may be supplied as pellets or granules. Do not touch the pellets or granules since sodium hydroxide is very caustic and can cause chemical burns. Use a scoop or forceps and wear gloves. Sodium hydroxide should be weighed as rapidly as possible since it is hygroscopic and draws moisture from the air, becoming sticky. Keep the lid on the bottle as much of the time as possible. If the sodium hydroxide is in pellet form, you will not be able to weigh out an exact amount. Don't try to break up a pellet. When making up the sodium hydroxide solution, you should always add sodium hydroxide to water. Adding water to sodium hydroxide can cause an intense localized temperature increase that can boil the liquid and cause splattering.

Prepare one of the three soaps given below. For each soap, you will need to calculate the amounts of sodium hydroxide, fat or oil, and water required. You should use the table given in Experiment 26A and the calculation examples given in the introduction to Experiment 26A as a guide. The calculations shown are based on a combined total of 100 g of fat or oil in the recipe. Reduce this by two-thirds for this experiment. The easiest way to do this is to first assume a total of 100 g of the fat/oil. Calculate the amount of base, fat or oil, and water as shown in the examples, then divide all of your values by 3. Note that in Soap #1 and #3, the superfat percentage is different than the calculation examples. You will need to make the appropriate adjustments. **Check your calculations with your instructor before you weigh out any of the ingredients.**

To prepare your soap, follow the instructions given in Experiment 26A with your calculated quantities for the recipe of your choice.

Besides using different fat and oil combinations, adding a scent is a great way to customize your soap. Since smell is closely tied to emotional memory, even a little bit goes a long way. To add a scent(s) to your soap, add 1–2 g of essential oil(s) or fragrance oil(s) of your choosing when you reach the “trace” stage or at the end of the 60-minute reaction period. Mix thoroughly for 5 minutes before pouring into the mold.

Before using your soap or taking your soap home, you must test the pH to make sure that it is not too basic. Refer to Experiment 26A for specific instructions on how to test the pH. Test the pH of the soap after it has set in the mold for about 1 hour, and again after 1 week, 2 weeks, and 3 weeks. If the pH is still high (>10) after 3 weeks, discard the soap since it could cause harm to your body. You should include the following table with your entries for the pH values in your report.

Time Interval	Initial pH	1 week	2 weeks	3 weeks
pH Value				

Soap 1: 100% Coconut Oil Soap, 20% Superfat

Description: Great one oil body bar that breaks the rules of soap making.

Comments: This soap should firm up within 1 hour after pouring into the mold.

Soap 2: 30% Coconut Oil, 70% Canola oil, 5% Superfat

Description: Basic all-around vegetarian soap.

Comments: This soap may require 3–4 hours to firm up after pouring into the mold.

Soap 3: 40% Vegetable Shortening (Crisco), 30% Coconut Oil, 30% Olive Oil, 8% Superfat

Description: Awesome vegetarian bath soap with added moisturizing quality.

Comments: This soap may require 12 hours to firm up after pouring into the mold.

QUESTIONS

1. Contrast the reaction to form biodiesel with the reaction to form soap. Why do they form different products given that both reactions start out with a base and a fat or oil?
2. Draw the structure of a fat or oil molecule that contains 1/3 lauric acid, 1/3 palmitic acid, and 1/3 oleic acid in its composition. Is the structure you drew the only unique answer possible? Explain.
3. Do a specific calculation of the amount of NaOH required to make a soap that is made by saponifying a mixture containing 20% avocado oil, 30% lard, and 50% canola oil. Assume that the total weight of the three fats/oils is 100 g and a superfat of 5%.
4. Give possible structures of the main components in a soap where the fat or oil was saponified with only two-thirds of the correct amount of sodium hydroxide. Use a generic fat or oil where the alkyl portions are represented by an R group.
5. Detergents are often used rather than soaps. Using an Internet search, find the structures for several detergents. What advantages do detergents have over soaps?

Petroleum and Fossil Fuels

Crude petroleum is a liquid that consists of hydrocarbons, as well as some related sulfur, oxygen, and nitrogen compounds. Other elements, including metals, may be present in trace amounts. Crude oil is formed by the decay of marine animal and plant organisms that lived millions of years ago. Over many millions of years, under the influence of temperature, pressure, catalysts, radioactivity, and bacteria, the decayed matter was converted into what we now know as crude oil. The crude oil is trapped in pools beneath the ground by various geological formations.

Most crude oils have a specific gravity between 0.78 and 1.00 g/mL. As a liquid, crude oil may be as thick and black as melted tar or as thin and colorless as water. Its characteristics depend on the particular oil field from which it comes. Pennsylvania crude oils are high in straight-chain alkane compounds (called **paraffins** in the petroleum industry); those crude oils are therefore useful in the manufacture of lubricating oils. Oil fields in California and Texas produce crude oil with a higher percentage of cycloalkanes (called **naphthenes** by the petroleum industry). Some Middle East fields produce crude oil containing up to 90% cyclic hydrocarbons. Petroleum contains molecules in which the number of carbons ranges from 1 to 60.

When petroleum is refined to convert it into a variety of usable products, it is initially subjected to a fractional distillation. Table 1 lists the various fractions obtained from fractional distillation. Each of these fractions has its own particular uses. Each fraction may be subjected to further purification, depending on the desired application.


The gasoline fraction obtained directly from the distillation of crude oil is called **straight-run gasoline**. An average barrel of crude oil will yield about 19% straight-run gasoline. This yield presents two immediate problems. First, there is not enough gasoline contained in crude oil to satisfy current needs for fuel to power automobile engines. Second, the straight-run gasoline obtained from crude oil is a poor fuel for modern engines. It must be “refined” at a chemical refinery.

The initial problem of the small quantity of gasoline available from crude oil can be solved by **cracking** and **polymerization**. Cracking is a refinery process by which large hydrocarbon molecules are broken down into smaller molecules. Heat


Table 1 Fractions Obtained from the Distillation of Crude Oil

<i>Petroleum Fraction</i>	<i>Composition</i>	<i>Commercial Use</i>
Natural gas	C ₁ to C ₄	Fuel for heating
Gasoline	C ₅ to C ₁₀	Motor fuel
Kerosene	C ₁₁ to C ₁₂	Jet fuel and heating
Light gas oil	C ₁₃ to C ₁₇	Furnaces, diesel engines
Heavy gas oil	C ₁₈ to C ₂₅	Motor oil, paraffin wax, petroleum jelly
Residuum	C ₂₆ to C ₆₀	Asphalt residual oils, waxes


and pressure are required for cracking, and a catalyst must be used. Silica–alumina and silica–magnesia are among the most effective cracking catalysts. A mixture of saturated and unsaturated hydrocarbons is produced in the cracking process. If gaseous hydrogen is also present during the cracking, only saturated hydrocarbons are produced. The hydrocarbon mixtures produced by these cracking processes tend to have a fairly high proportion of branched-chain isomers. These branched isomers improve the quality of the fuel.

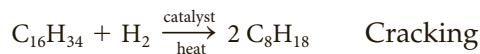


In the polymerization process, also carried out at a refinery, small molecules of alkenes are caused to react with one another to form larger molecules, which are also alkenes.

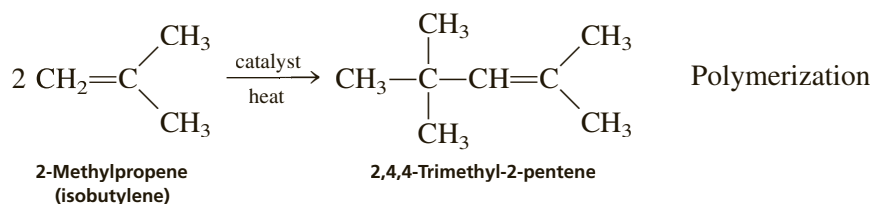

The newly formed alkenes may be hydrogenated to form alkanes. The reaction sequence shown here is a very common and important one in petroleum refining because the product, 2,2,4-trimethylpentane (or “isooctane”), forms the basis for determining the quality of gasoline. By these refining methods, the percentage of gasoline that can be obtained from a barrel of crude oil may rise to as much as 45% or 50%.

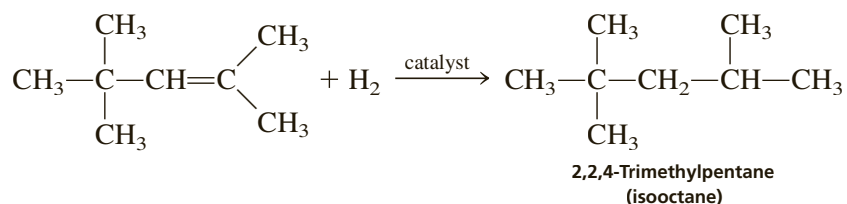

The internal combustion engine, as it is found in most automobiles, operates in four cycles or **strokes**. They are illustrated in the figure. The power stroke is of

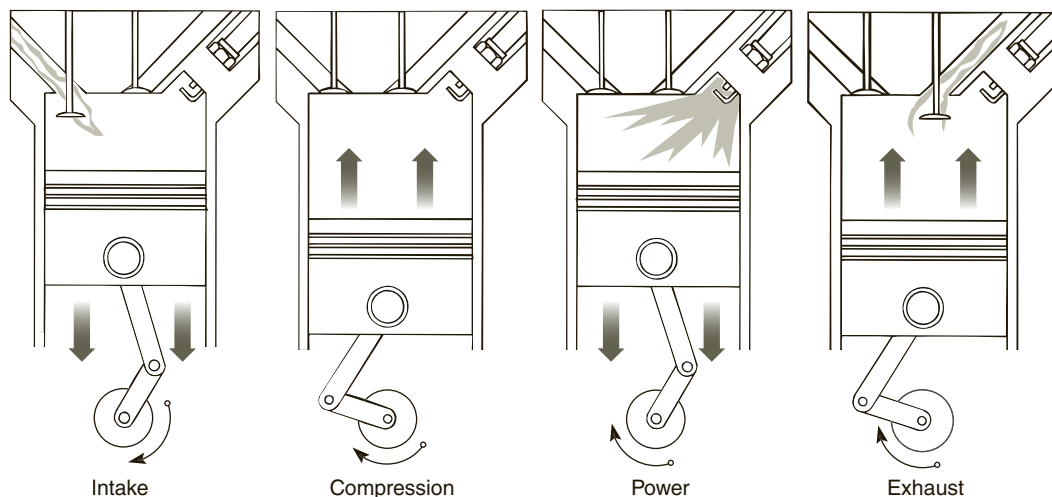

Operation of a four-cycle engine.

greatest interest from the chemical point of view because combustion occurs during this stroke.

When the air–fuel mixture is ignited, it does not explode. Rather, it burns at a controlled, uniform rate. The gases closest to the spark are ignited first; then they in turn ignite the molecules farther from the spark; and so on. The combustion proceeds in a wave of flame or a **flame front**, which starts at the spark plug and proceeds uniformly outward from that point until all the gases in the cylinder have been ignited. Because a certain time is required for this burning, the initial spark is timed to ignite just before the piston has reached the top of its travel. In this way, the piston will be at the very top of its travel at the precise instant that the flame front and the increased pressure that accompanies it reach the piston. The result is a smoothly applied force to the piston, driving it downward.

If heat and compression should cause some of the air–fuel mixture to ignite before the flame front has reached it or to burn faster than expected, the timing of the combustion sequence is disturbed. The flame front arrives at the piston before the piston has reached the very top of its travel. When the combustion is not perfectly coordinated with the motion of the piston, we observe **knocking**, or **detonation** (sometimes called “pinging”). The transfer of power to the piston under these conditions is much less effective than in normal combustion. The wasted energy is merely transferred to the engine block in the form of additional heat. The opposing forces that occur in knocking may eventually damage the engine.




The tendency of a fuel to knock is a function of the structures of the molecules composing the fuel. Normal hydrocarbons, those with straight carbon chains, have a greater tendency to lead to knocking than do alkanes with highly branched chains. A fuel can be classified according to its antiknock characteristics. The most important rating system is the octane rating of gasoline. In this method of classification, the antiknock properties of a fuel are compared in a test engine with the antiknock properties of a standard mixture of heptane and 2,2,4-trimethylpentane. This latter compound is called “isooctane,” hence the name *octane rating*. A fuel that has the same antiknock properties as a given mixture of heptane and isooctane has an octane rating numerically equal to the percentage of isooctane in that reference mixture. Today’s 87-octane unleaded gasoline is a mixture of compounds that have, taken together, the same antiknock characteristics as a test fuel composed of 13% heptane and 87% isooctane. Other substances besides hydrocarbons may also have high resistance to knocking. Table 2 presents a list of organic compounds with their octane ratings.

Table 2 Octane Ratings of Organic Compounds

<i>Compound</i>	<i>Octane Number</i>	<i>Compound</i>	<i>Octane Number</i>
Octane	–19	1-Butene	97
Heptane	0	2,2,4-Trimethylpentane	100
Hexane	25	Cyclopentane	101
Pentane	62	Ethanol	105
Cyclohexane	83	Benzene	106
1-Pentene	91	Methanol	106
2-Hexene	93	Methyl <i>tert</i> -butyl ether	116
Butane	94	<i>m</i> -Xylene	118
Propane	97	Toluene	120


Note: The octane values in this table are determined by the **research method**.

Several chemical refining processes are used to improve the octane rating of gasoline and to increase the percentage of gasoline that can be obtained from petroleum. Some of these reactions, collectively known as **reforming**, are dehydrogenation, dealkylation, cyclization, and isomerization. The products of these reactions, sometimes referred to as **reformates**, contain many branched alkanes and aromatic compounds. Several examples of reforming reactions are:

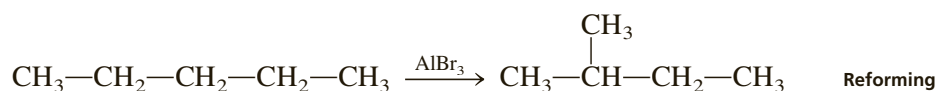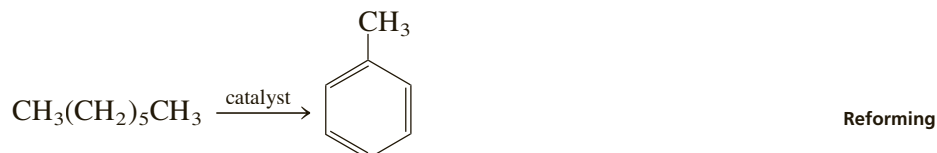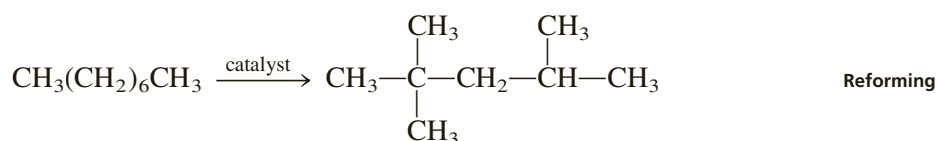


Other chemical reactions, referred to as **alkylation**, can also be used to increase the octane rating. Alkylation involves the catalytic addition of an alkane to an alkene, such as 2-methylpropane to propene or butane. The products of these reactions are sometimes referred to as **alkylates**. Another refining process, called **hydrocracking** (cracking in the presence of hydrogen gas), also produces hydrocarbons that reduce knocking.

None of these processes converts all the normal hydrocarbons into branched-chain isomers; consequently, additives are also put into gasoline to improve the octane rating of the fuel. Before 1996, the most common additive used to reduce knocking has been **tetraethyllead**. Gasoline that contains tetraethyllead is called **leaded** gasoline, whereas gasoline produced without tetraethyllead is sometimes called **unleaded** gasoline. Because of concern over the possible health hazard associated with emission of lead into the atmosphere the Environmental Protection Agency began in 1973 to limit the amount of tetraethyllead in gasoline. In 1996, the Clean Air Act completely banned the sale of leaded gasoline for use in all on-road vehicles. Many other countries have followed with similar bans; however, some countries in Eastern Europe, the Middle East, and Africa continue to use leaded gasoline.

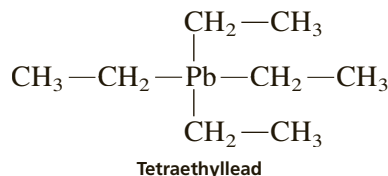


To replace tetraethyllead, oil companies have developed other additives and strategies that will improve the octane rating of gasoline without producing harmful emissions. One approach is to increase quantities of hydrocarbons that have very high antiknock properties themselves. Typical are the aromatic hydrocarbons, including benzene, toluene, and xylene. Such compounds are natural components of most crude petroleum, and additional aromatic compounds can be added to gasoline to improve the quality. Increasing the proportion of aromatic hydrocarbons

brings with it certain hazards, however. These substances are toxic, and benzene is considered a serious carcinogenic hazard. The risk that illness will be contracted by workers in refineries, and especially by persons who work in service stations, is increased. A safer approach is to increase the amount of alkylates.

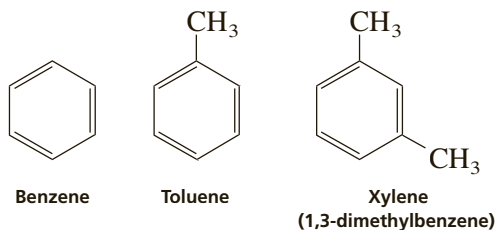

Research has also been directed toward development of nonhydrocarbon compounds that can improve the quality of unleaded gasoline. To this end, compounds such as methyl *tert*-butyl ether (MTBE), ethanol, and other **oxygenates** (oxygen-containing compounds) are added to improve the octane rating of fuels. Ethanol is attractive because it is formed by fermentation of living material, a renewable resource (see essays “Biofuels,” which precedes Experiment 28, and “Ethanol and Fermentation Chemistry,” which precedes Experiment 18). Ethanol not only would improve the antiknock properties of gasolines, but also would potentially help the country to reduce its dependence on imported petroleum. Substituting ethanol for hydrocarbons in petroleum would have the effect of increasing the “yield” of fuel produced from a barrel of crude oil. As in many stories that are too good to be true, it is not clear that the energy needed to produce the ethanol by fermentation and distillation is significantly smaller than the amount of energy that is produced when the ethanol is burned in an engine!

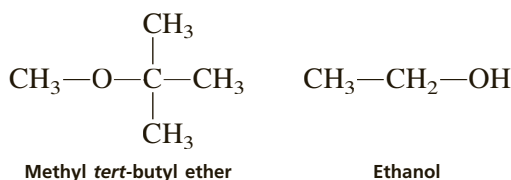


In an effort to improve air quality in urban areas, the Clean Air Act of 1990 mandated the addition of oxygen-containing compounds in many urban areas during the winter (November to February). These compounds are expected to reduce carbon monoxide emissions produced when the gasoline burns in cold engines by helping to oxidize carbon monoxide to carbon dioxide. They also help to reduce the amount of ozone created by emission products reacting in sunlight; and they increase the octane rating. Refineries add “oxygenates,” such as ethanol or methyl *tert*-butyl ether, to the gasoline sold in the carbon monoxide-containment areas. By law, gasoline must contain at least 2.7% oxygen by weight, and the areas must use it for a minimum of the four winter months. In 1995, the Clean Air Act also required that **reformulated gasoline** (RFG) be sold year round in sites with the worst ground-level ozone concentrations. RFG must contain a minimum of 2% oxygen by weight.

Although methyl *tert*-butyl ether is still used in some states, the use of ethanol is much more common. There are several reasons for the preference for ethanol. First, ethanol is cheaper than MTBE because of special tax breaks and subsidies that have been granted to producers of ethanol formed by fermentation. Second, there has been much concern that MTBE may cause health problems, and there have been some widely publicized occurrences of groundwater contamination by

MTBE. Furthermore, people notice the odor of gasoline more easily when MTBE is present in the fuel. Because of these concerns, the use of MTBE was outlawed by California in January 2004, and many other states have issued similar or partial bans. It is possible that a complete ban on MTBE in the United States will follow. Therefore, ethanol has become the preferred oxygenate for gasoline. However, there are disadvantages with the use of ethanol, too. There is some evidence that because ethanol is more volatile than MTBE, it may increase the emission of chemicals such as volatile organic compounds (VOCs) that contribute to smog. This is a concern especially during the warmer months. In addition, studies have suggested that fuel with ethanol increases the formation of atmospheric acetaldehyde. Because acetaldehyde is a precursor to peroxyacetyl nitrate, it is possible that *increased* air pollution results from use of ethanol as an oxygenate. Other oxygenates such as ethyl *tert*-butyl ether and methanol are also being considered.

The number of grams of air required for the complete combustion of one mole of gasoline (assuming the formula C_8H_{18}) is 1.735 grams. This gives rise to a theoretical air–fuel ratio of 15.1:1 for complete combustion. For several reasons, however, it is neither easy nor advisable to supply each cylinder with a theoretically correct air–fuel mixture. The power and performance of an engine improve with a slightly richer mixture (lower air–fuel ratio). Maximum power is obtained from an engine when the air–fuel ratio is near 12.5:1, and maximum economy is obtained when the air–fuel ratio is near 16:1. Under conditions of idling or full load (that is, acceleration), the air–fuel ratio is lower than what would be theoretically correct. As a result, complete combustion does not take place in an internal combustion engine, and carbon monoxide (CO) is produced in the exhaust gases. Other types of nonideal combustion behavior give rise to the presence of unburned hydrocarbons in the exhaust. The high combustion temperatures cause the nitrogen and oxygen of the air to react, forming a variety of nitrogen oxides in the exhaust. Each of these materials contributes to air pollution. Under the influence of sunlight, which has enough energy to break covalent bonds, these materials may react with each other and with air to produce **smog**, which contributes to many health problems. Smog consists of **ozone**, which deteriorates rubber and damages plant life; **particulate matter**, which produces haze; **oxides of nitrogen**, which produce a brownish color in the atmosphere; and a variety of eye irritants, such as **peroxyacetyl nitrate** (PAN). Sulfur compounds in the gasoline may lead to the production of noxious sulfur-containing gases in the exhaust.

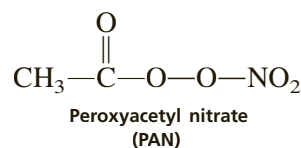

Efforts to reverse the trend of deteriorating air quality caused by automotive exhaust have taken many forms. The advent of catalytic converters, which are mufflerlike devices containing catalysts that can convert carbon monoxide, unburned hydrocarbons, and nitrogen oxides into harmless gases, has resulted from such efforts. Some success in reducing exhaust emissions has been attained by modifying the design of combustion chambers of internal combustion engines. Additionally, the use of computerized control of ignition systems has achieved positive results.

It should be obvious from this discussion that there are many factors considered in the formulation of gasoline. The gasoline produced today consists of several hundred compounds! There is substantial variation in the actual composition,

depending on the local climate and regional environmental regulations. The approximate composition is 15% C₄–C₈ straight-chain alkanes, 25% C₄–C₁₀ branched alkanes, 10% cycloalkanes, less than 25% aromatic compounds, and 10% straight-chain and cyclic alkenes.

Although much has been accomplished in terms of making it safer to use gasoline, there is another looming problem having to do with the supply of petroleum. The amount of petroleum and other fossil fuels in the world is finite. In 1956, Marion King Hubbert, a Shell Oil geophysicist, predicted that the U.S. production of oil would reach a peak around 1970, and from then on the amount extracted would decline significantly. Although most people ignored his warning, it did peak at 9 million barrels per day in 1970 and has been declining ever since, with about 6 million barrels per day being produced in 2004. Many experts have used similar methods of analysis to make predictions about when the world's supply of oil will peak; and although there is much variation in the actual year predicted, most experts agree that the peak has already occurred. Because the demand for petroleum continues to increase every year, it is clear that declining petroleum production will have a dramatic effect on how we live. Not only is petroleum the main source of fuel used for transportation but it also provides the raw materials for a wide variety of other products, including plastics, drugs, and pesticides. Although it is possible that the decrease in production of petroleum may be partially offset by more dependence on natural gas and coal, the amount of these fossil fuels is also finite, and it seems inevitable that major adjustments will need to be made as the availability of fossil fuels declines.

Many developments in recent years have addressed some of the emission problems associated with burning gasoline and the need to stretch the supply of fossil fuels. These developments involve changes in the design of automobile engines and in the use of different fuels.

Some of the success in reducing exhaust emission has been attained by modifying the design of combustion chambers of internal combustion engines. Additionally, the use of computerized control of ignition systems has helped to reduce the level of pollutants emitted. Another strategy that could be implemented without any technological changes would be to increase fuel standard requirements, thus improving the average miles per gallon. Because this would result in less gasoline consumption, there would also be less emission of pollutants.

Diesel engines have been used in automobiles for more than 20 years. These engines require a different fraction of crude oil (see Table 1 at the beginning of this essay) than gasoline, and they have been improved significantly since the initial highly polluting diesel vehicles. The diesel engine has the advantage of producing only small quantities of carbon monoxide and unburned hydrocarbons. It does, however, produce significant amounts of nitrogen oxides, soot (containing polynuclear aromatic hydrocarbons), and odor-causing compounds. Presently, the emission standards for diesel automobiles are more lenient than for those burning gasoline. More stringent standards were implemented in 2006 and 2009. Diesel automobiles yield higher fuel mileage than gasoline engines of a similar size; however, more oil must be refined in order to produce diesel fuel compared to gasoline. In the United States, about 3% of all new automobiles have diesel engines, whereas in Europe, about 40% of the new automobiles sold are diesel. Biodiesel, which is a chemically altered vegetable oil that can even be produced in one's garage using discarded cooking oil, can also be used in today's diesel engines and results in fewer harmful pollutants compared to regular diesel fuel. However, the mileage is slightly less, and it would not be possible to produce enough of this fuel for more than a small percentage of the cars on the road today.

Another possible fuel is methanol, which is produced from natural gas, coal, or biomass. Studies indicate that the amount of principal pollutants in automobiles is lowered when methanol is used instead of gasoline, but methanol is more corrosive and extensive engine modifications must be made. Other fuels that show promise are hydrogen, methane (natural gas), and propane; however, storage and delivery of these fuels, which are gases at room temperature, are more difficult and other significant technical problems also must be solved.

It is now clear that the most significant problem related to the combustion of fossil fuels is likely global warming, due to the increasing concentration of carbon dioxide in the atmosphere. Most of the radiant energy from the sun passes through the earth's atmosphere and reaches the earth, where much of this energy is converted into heat. Most of this heat in the form of infrared radiation is radiated away from the earth. Carbon dioxide and other atmospheric compounds, such as water and methane, can absorb this infrared radiation. When this heat energy is released by these molecules, it radiates in all directions—including back toward the earth. The retention of some of this heat is referred to as the **greenhouse effect**. The greenhouse effect is extremely valuable in terms of keeping the temperature of the earth in a range where life can exist. However, the temperature of the earth has been increasing during the past century, likely because of the increase in the amount of carbon dioxide in the atmosphere. Most of this additional carbon dioxide is produced by the combustion of fossil fuels. There is much concern that if the temperature of the earth continues to increase, the implication for life on the earth could be devastating. The sea levels might rise high enough to force millions of people living in coastal areas to migrate, and the negative effect on farming and fresh water sources could have a serious impact on people in all parts of the world.

Hybrid-electric automobiles have become an attractive alternative to the standard automobile in the United States. Hybrid cars combine a small fuel-efficient combustion engine with an electric motor and battery. The electric motor can assist the gas engine when more power is needed, and the battery is recharged while the car is slowing down or coasting. This results in greater fuel efficiency, as well as a drastic reduction in the amount of carbon dioxide released and smog-forming pollutants. Even greater fuel efficiency is possible with diesel hybrid cars that are now being developed. In the past few years, there has been increasing interest in the development of electric plug-in cars that run off a large storage battery. These batteries can be charged at night when the overall electrical demand on the grid is low, and the cars can be driven about 30–200 miles on a charge, depending on the type of battery. If the electricity were generated by a renewable energy source such as solar, wind, or geothermal, then the contribution to the greenhouse effect by driving electric cars would be minimal.

Another recent promising development is the use of fuel cells that can produce electrical energy from hydrogen. This electrical energy is then used by an electric motor to propel the automobile. Although there are many proponents of hydrogen fuel cells who believe that this technology can play a major role in reducing our dependency on fossil fuels, significant technological challenges must first be overcome. The task of developing a hydrogen energy infrastructure would also be costly. Furthermore, most of the hydrogen now produced comes from natural gas or coal, and this process also requires energy.

It should be clear that the use of fossil fuels poses many challenges and opportunities. How we utilize fossil fuels will change in the next few decades, and chemistry will play a significant role.

REFERENCES

- Ashley, S. On the Road to Fuel-Cell Cars. *Sci. Am.* **2005**, 292 (Mar), 62.
- Chen, C. T. Understanding the Fate of Petroleum Hydrocarbons in the Subsurface Environment. *J. Chem. Educ.* **1992**, 69 (May), 357.
- Goodstein, D. *Out of Gas: The End of the Age of Oil*; W. W. Norton & Company Inc.: New York, 2004.
- Hogue, C. Rethinking Ethanol. *Chem. Eng. News* **2003**, 81 (Jul 28), 28.
- Illman, D. Oxygenated Fuel Cost May Outweigh Effectiveness. *Chem. Eng. News* **1993**, 71 (Apr 12), 28.
- Jacoby, M. Fuel Cells Move Closer to Market. *Chem. Eng. News* **2003**, 81 (Jan 20), 328.
- Kimmel, H. S.; Tomkins, R. P. T. A Course on Synthetic Fuels. *J. Chem. Educ.* **1985**, 62 (Mar), 249.
- Mielke, H. W. Lead in the Inner Cities. *Am. Sci.* **1999**, 87 (Jan), 63–73.
- Monahan, P.; Friedman, D. Diesel or Gasoline? Fuel for Thought. *Catalyst* **2004**, 3 (Spring), 1.
- Ritter, S. Gasoline. *Chem. Eng. News* **2005**, 83 (Feb 21), 37.
- Schriesheim, A.; Kirschenbaum, I. The Chemistry and Technology of Synthetic Fuels. *Am. Sci.* **1981**, 69 (Sep–Oct), 536.
- Seyferth, D. The Rise and Fall of Tetraethyllead. 1. Discovery and Slow Development in European Universities, 1853–1920. *Organometallics* **2003**, 22 (Jun 9), 2346–2357.
- Seyferth, D. The Rise and Fall of Tetraethyllead. 2. *Organometallics* **2003**, 22 (Dec 8), 5154–5178.
- Shreve, R. N.; Brink, J. Petrochemicals. *The Chemical Process Industries*, 4th ed.; McGraw-Hill: New York, 1977.
- Shreve, R. N.; Brink, J. Petroleum Refining. *The Chemical Process Industries*, 4th ed.; McGraw-Hill: New York, 1977.
- U.S. Environmental Protection Agency. EPA Takes Final Step in Phaseout of Leaded Gasoline. <http://www.epa.gov/history/topics/lead/0.2.htm> (accessed Jan 29, 1996).
- U.S. Environmental Protection Agency. MTBE in Fuels. <http://www.epa.gov/mtbe/gas.htm> (accessed Jun 26, 2005).
- Vartanian, P. F. The Chemistry of Modern Petroleum Product Additives. *J. Chem. Educ.* **1991**, 68 (Dec), 1015.
- Wald, M. L. Questions About a Hydrogen Economy. *Sci. Am.* **2004**, 291 (May), 62.

Gas-Chromatographic Analysis of Gasolines

Gasoline

Gas chromatography

In this experiment, you will analyze samples of gasoline by gas chromatography. From your analysis, you should learn something about the composition of these fuels. Although all gasolines are compounded from the same basic hydrocarbon components, each company blends these components in different proportions to obtain a gasoline with properties similar to those of competing brands.

Sometimes the composition of the gasoline may vary, depending on the composition of the crude petroleum from which the gasoline was derived. Frequently, refineries vary the composition of gasoline in response to differences in climate, seasonal changes, or environmental concerns. In the winter or in cold climates, the relative proportion of butane and pentane isomers is increased to improve the volatility of the fuel. This increased volatility permits easier starting. In the summer or in warm climates, the relative proportion of these volatile hydrocarbons is reduced. The decreased volatility reduces the possibility of forming a vapor lock. Occasionally, differences in composition can be detected by examining the gas chromatograms of a particular gasoline over several months. In this experiment, we will not try to detect such small differences.

There are different octane rating requirements for “regular” and “premium” gasolines. You may be able to observe differences in the composition of these two types of fuels. You should pay particular attention to increases in the proportions of those hydrocarbons that raise octane ratings in the premium fuels.

In some areas of the country, manufacturers are required from November to February to control the amounts of carbon monoxide produced when the gasoline burns. To do this, they add oxygenates, such as ethanol or methyl *tert*-butyl ether (MTBE), to the gasoline. You should try to observe the presence of these oxygenates, which may be observed in gasolines produced in carbon monoxide–containment areas. Because MTBE has been banned or partially banned in most states (see previous essay), it is unlikely that you will observe MTBE.

The class will analyze samples of the unleaded regular and premium gasolines in use today. If available, the class will analyze oxygenated fuels. If different brands are analyzed, equivalent grades from the different companies should be compared.

Discount service stations usually buy their gasoline from one of the large petroleum-refining companies. If you analyze gasoline from a discount service station, you may find it interesting to compare that gasoline with an equivalent grade from a major supplier, noting particularly the similarities.

REQUIRED READING

New: Technique 22 Gas Chromatography
Essay Petroleum and Fossil Fuels

SPECIAL INSTRUCTIONS

Your instructor may want each student in the class to obtain a sample of gasoline from a service station. The instructor will compile a list of the different gasoline companies represented in the nearby area. Each student will then be assigned to collect a sample from a different company. You should collect the gasoline sample in a labeled screw-cap jar. An easy way to collect a gasoline sample for this experiment is to drain the excess gasoline from the nozzle and hose into the jar after the gasoline tank of a car has been filled. The collection of gasoline in this manner must be done *immediately after* the gas pump has been used. If not, the volatile components of the gasoline may evaporate, thus changing the composition of the gasoline. Only a small sample (a few *milliliters*) is required because the gas-chromatographic analysis requires no more than a few *microliters* (μL) of material. Be certain to close the cap of the sample jar tightly to prevent the selective evaporation of the most volatile components. The label on the jar should list the brand of gasoline and the grade (regular, premium, oxygenated, etc.). Alternatively, your instructor may supply samples for you.

CAUTION



Gasoline contains many highly volatile and flammable components. Do not breathe the vapors, and do not use open flames near gasoline.

This experiment may be assigned along with another short one because it requires only a few minutes of each student's time to carry out the actual gas chromatography. For this experiment to be conducted as efficiently as possible, you may be asked to schedule an appointment for using the gas chromatograph.

SUGGESTED WASTE DISPOSAL

Dispose of all gasoline samples in the container designated for nonhalogenated wastes.

NOTES TO THE INSTRUCTOR

You need to adjust your gas chromatograph to the proper conditions for the analysis. We recommend that you prepare and analyze the reference mixture listed in the Procedure section. Most chromatographs will be able to separate this mixture cleanly with the possible exception of the xylenes. One possible set of conditions for a Gow-Mac model 69-350 chromatograph is the following: column temperature, 110–115°C; injection port temperature, 110–115°C; carrier gas flow rate, 40–50 mL/min; column length, approximately 12 ft. The column should be packed with a nonpolar stationary phase similar to silicone oil (SE-30) on Chromosorb W or with some other stationary phase that separates components principally according to boiling point.

The chromatograms shown in this experiment were obtained on a Hewlett Packard model 5890 gas chromatograph. A 30-meter, DB 5 capillary column (0.32 mm, with 0.25 micron film) was used. A temperature program was run starting at 5°C and ramping to 105°C. Each run took about 8 minutes. A flame-ionization detector was used. The conditions are given in the Instructor's Manual. Superior separations are obtained using capillary columns, which are recommended. Even better results are obtained with longer columns.

PROCEDURE

Reference Mixture

First, analyze a standard mixture that includes pentane, hexane (or hexanes), benzene, heptane, toluene, and xylenes (a mixture of *meta*, *para*, and *ortho* isomers). Inject a 0.5- μ L sample or an alternative sample size as indicated by your instructor into the gas chromatograph. Measure the retention time of each component in the reference mixture on your chromatogram (see Technique 22, Section 22.7). The previously listed compounds elute in the order given (pentane first and xylenes last). Compare your chromatogram to the one posted near the gas chromatograph or the one reproduced in this experiment.

Your instructor or a laboratory assistant may prefer to perform the sample injections. The special microliter syringes used in the experiment are delicate and expensive. If you are performing the injections yourself, be sure to obtain instruction beforehand.

Oxygenated Fuel Reference Mixture

Oxygenated compounds are added to gasolines in carbon monoxide-containment areas during November through February. Currently, ethanol is used most commonly. It is much less likely that methyl *tert*-butyl ether will be found. Your instructor may have available a reference mixture that includes all the previously listed compounds and either ethanol or methyl *tert*-butyl ether. Again, you need to inject a sample of this mixture and analyze the chromatogram to obtain the retention times for each component in this mixture.

Gasoline Samples

Inject a sample of a regular, premium, or oxygenated gasoline into the gas chromatograph and wait for the gas chromatogram to be recorded. Compare the chromatogram to that of the reference mixture. Determine the retention times for the major components and identify as many of the components as possible. For comparison, gas chromatograms of a premium gasoline and the reference mixture are shown on the following page. On the list of the major components in gasolines, notice that the oxygenate methyl *tert*-butyl ether appears in the C₆ region. Does your oxygenated fuel show this component? See if you can notice a difference between regular and premium gasolines.

Analysis

Be certain to compare carefully the retention times of the components in each fuel sample with the standards in the reference mixture. Retention times of compounds vary with the conditions under which they are determined. It is best to analyze the reference mixture and each of the gasoline samples in succession to reduce the variations in retention times that may occur over time. Compare the gas chromatograms with those of students who have analyzed gasolines from other dealers.


Report

The report to the instructor should include the actual gas chromatograms, as well as an identification of as many of the components in each chromatogram as possible.


Major components in gasolines*

C ₄ Compounds	Isobutane
	Butane
C ₅ Compounds	Isopentane
	Pentane
	2,3-Dimethylbutane
C ₆ Compounds and oxygenates	2-Methylpentane
	3-Methylpentane
	Hexane
	Methyl <i>tert</i> -butyl ether (oxygenate)
	2,4-Dimethylpentane
C ₇ Compounds and aromatics (benzene)	Benzene (C ₆ H ₆)
	2-Methylhexane
	3-Methylhexane
	Heptane
	2,2,4-Trimethylpentane (isooctane)
	2,5-Dimethylhexane
C ₈ Compounds and aromatics (toluene, ethylbenzene, and xylenes)	2,4-Dimethylhexane
	2,3,4-Trimethylpentane
	2,3-Dimethylhexane
	Toluene (C ₇ H ₈)
	Ethylbenzene (C ₈ H ₁₀)
	<i>m</i> -, <i>p</i> -, <i>o</i> -Xylenes (C ₈ H ₁₀)
	1-Ethyl-3-methylbenzene
C ₉ Aromatic compounds	1,3,5-Trimethylbenzene
	1,2,4-Trimethylbenzene
	1,2,3-Trimethylbenzene

*Approximate order of elution.

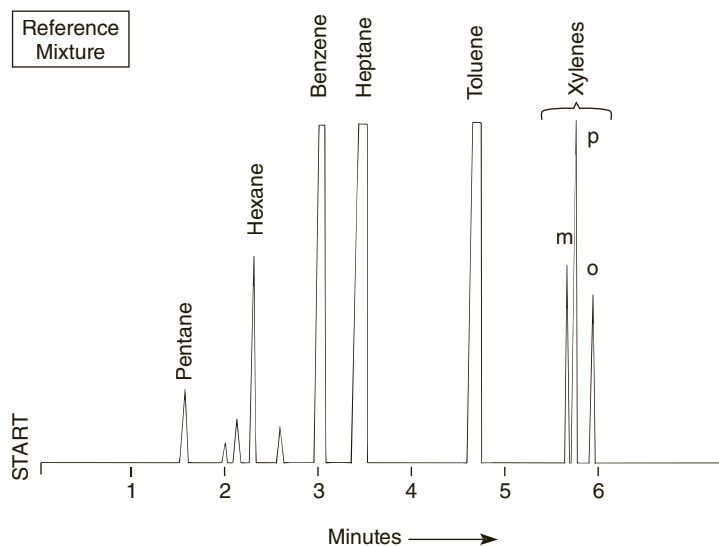

Gas chromatogram of the reference mixture.

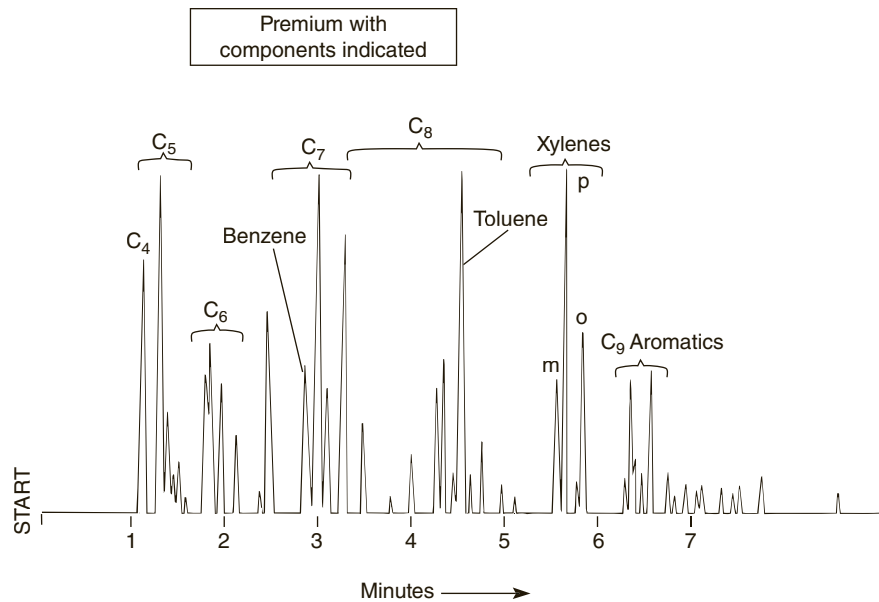

Gas chromatogram of a premium gasoline.

QUESTIONS

1. If you had a mixture of benzene, toluene, and *m*-xylene, what would be the expected order of retention times? Explain.
2. If you were a forensic chemist working for the police department and the fire marshal brought you a sample of gasoline found at the scene of an arson attempt, could you identify the service station at which the arsonist purchased the gasoline? Explain.
3. How could you use infrared spectroscopy to detect the presence of ethanol in an oxygenated fuel?

Biofuels


In recent years there has been an increasing interest in **biofuels**, fuels that are produced from biological materials such as corn or vegetable oil. These sources of biofuels are considered to be renewable because they can be produced in a relatively short time. On the other hand, **fossil fuels** are formed by the slow decay of marine animal and plant organisms that lived millions of years ago. Fossil fuels, which include petroleum, natural gas, and coal, are considered to be nonrenewable.

The increased emphasis on biofuels is due primarily to the increasing cost and demand for liquid fuels such as gasoline and diesel, and our desire to be less dependent on foreign oil. In addition to increased demand, the higher cost of petroleum may be related to the peak oil theory, discussed in the essay on petroleum and fossil fuels. According to this theory, the amount of petroleum in the earth is finite; and at some point, the total amount of petroleum produced each year will begin to decrease. Many experts believe that we have either already reached the peak in oil production, or we will reach it within a few years.

In addition to biofuels, the use of many other types of alternative energy sources has been increasing in recent years. Alternative energy sources such as solar, wind, and geothermal are used primarily to produce electricity, and they cannot replace liquid fuels such as gasoline and diesel. As long as we continue to depend on automobiles and other vehicles with the current engine technology, we will need to produce more liquid fuels. Because of this, the demand to produce more biofuels is very great. In this essay, we will focus on the biofuels ethanol and biodiesel.

Ethanol

The knowledge of how to produce ethanol from grains has been around for many centuries (see the essay “Ethanol and Fermentation Chemistry” that precedes Experiment 18). Until recently, most of the ethanol produced by fermentation was used mainly in alcoholic beverages. In 1978, Congress passed the National Energy Act, which encouraged the use of fuels such as Gasohol, a blend of gasoline with at least 10% ethanol produced from renewable resources. Ethanol can be produced by the fermentation of sugars such as sucrose, which is found in sugar cane or beets. In this country, it is more common to use corn kernels as the feedstock to produce ethanol. Corn contains starch, a polymer of glucose that must first be broken down into glucose units. This is usually accomplished by adding a mixture of enzymes that catalyze the hydrolysis of starch into glucose. Other enzymes are then added to promote the fermentation of glucose into ethanol:

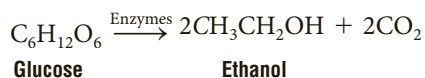

After fermentation, fractional distillation is used to separate the ethanol from the fermentation mixture. In Experiment 18, ethanol is separated from a fermentation mixture by fractional distillation.

The use of corn to produce ethanol as a biofuel has been strongly encouraged in the United States. Government subsidies have resulted in a higher production of corn in the Midwest, and many new ethanol refineries have also been built. However, it is now clear that use of ethanol as a biofuel has some significant drawbacks. First, as more corn is planted and used for fuel production, less corn and other crops are available as a source of food. This has led to food shortages and higher prices, which is especially hard on people who are already struggling to get enough food. Second, it now appears that the total amount of energy expended to grow corn and to produce ethanol is almost as much as the amount of energy released by burning the ethanol. Third, recent studies have indicated that growing corn to produce ethanol for use as a fuel results in the production of more greenhouse gases than the use of similar amounts of fossil fuels. Therefore, the use of corn ethanol may actually increase global warming compared to fossil fuels. In spite of these drawbacks, given that so much investment in corn ethanol has already been made, it is still likely that corn will continue to be a source of ethanol in this country for some time to come.


One alternative to corn ethanol is **cellulosic ethanol**. Sources of cellulose that can be used to produce ethanol include fast-growing grasses such as switchgrass, agricultural waste such as corn stalks, and waste wood from the milling of lumber. Like starch, cellulose is a polymer of glucose, but the structure is slightly different than starch and it is much more difficult to break down. Cellulose can be broken down by acid or base treatment at high temperature and by hydrolysis reactions with enzymes. Once the cellulose is broken down into glucose, it can be fermented to produce ethanol, just like with corn starch. Cellulosic ethanol addresses some of the drawbacks for corn ethanol mentioned in the previous paragraph. Many of the sources of cellulosic ethanol can be grown on non-arable land that would not normally be used to produce food. It also appears that the overall energy production is more favorable than with corn ethanol. Finally, the contribution to greenhouse gases is not so great. The first new cellulosic ethanol facilities began production in 2014.

Evaluating biofuels in terms of contribution to global warming is difficult to do. Initially, it was believed that all biofuels produced less greenhouse gases than fossil fuels. This is because carbon dioxide is absorbed by the plants as they grow, which helps to offset the carbon dioxide that is released when the biofuel is burned. However, recent studies suggest that the situation is more complicated. In order to grow the crops required to make biofuels and to replace the food crops that are now used to make biofuels, it is often necessary to destroy forestland. Forests are much more effective than farmland at absorbing carbon dioxide from the air.

Another option for producing ethanol exists that may have advantages over both of the methods described above. This newer option involves the conversion of carbon-containing matter into **syngas**. Almost any material that contains carbon, such as municipal waste, old tires, or agricultural waste, can be used. The feedstock is gasified into a mixture of carbon monoxide and hydrogen, which is known as syngas. Syngas can then be catalytically converted into ethanol. This process is much more efficient energetically than the methods described above and also creates less greenhouse gases, especially when the feedstock is some kind of waste material. Furthermore, these feedstocks do not compete with food crops.

Biodiesel

Another biofuel that is widely used in the United States is **biodiesel**. Biodiesel is produced from fats or oils in a base-catalyzed transesterification reaction:

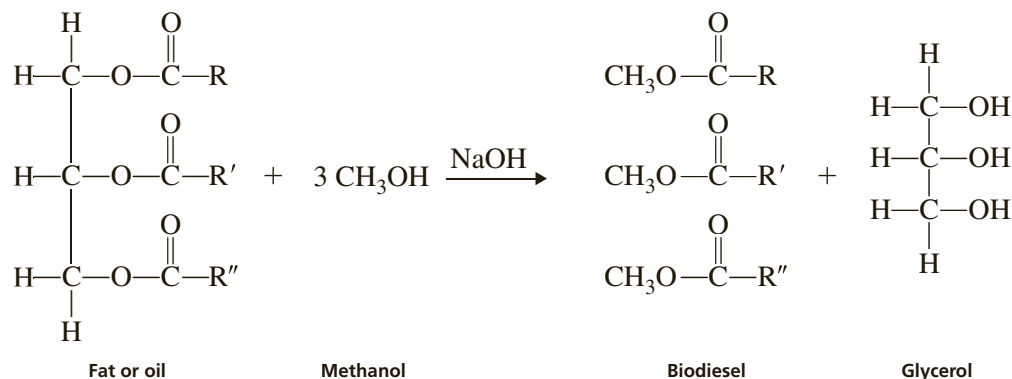

Because the R groups may have different numbers of carbons and double bonds, biodiesel is a mixture of different molecules, all of which are methyl esters of fatty acids. Most of the R groups have 12–18 carbons arranged in straight chains. Any kind of vegetable oil can be used to make biodiesel, but the most common ones used are the oils from soybean, canola, and palm. In Experiment 28, biodiesel is made from coconut oil and other vegetable oils.

Biodiesel has similar properties to the diesel fuel that is produced from petroleum, and it can be burned in any vehicle with a diesel engine or in furnaces that burn diesel fuel. It should be noted that vegetable oil can also be burned as a fuel, but because the viscosity of vegetable oil is somewhat greater than diesel fuel, engines must be modified in order to burn vegetable oil.

How does biodiesel compare with ethanol? Like corn ethanol, growing the vegetables required to produce the oil feedstock results in diverting farmland from growing food to producing fuels. In fact, this is more of a problem with biodiesel because more land is required to produce an equivalent amount of fuel compared to corn ethanol. The net energy produced by biodiesel is greater than for corn ethanol, but less than for cellulosic ethanol. Finally, it appears that the production of biodiesel, like ethanol, produces more greenhouse gases than fossil fuels, again because forested land must be destroyed in order to grow the vegetables required to produce biodiesel.

Some alternative approaches for making biodiesel exist that could address some of these issues. Algae can produce oils that can be used to make biodiesel. Algae can be grown in ponds or even waste water and do not require the use of farmland. The algae oil can be converted into biodiesel in the same way that vegetable oil is converted. Recently, a different chemical method for making biodiesel from vegetable oil has been developed. This method utilizes a sulfated zirconia catalyst that is placed in a column, similar to column chromatography. As the mixture of oil and alcohol is passed through the column at high temperature and pressure, biodiesel is produced and elutes from the bottom of the column. The process is much more efficient than the current methods used to produce biodiesel. An interesting side story related to this process is that the original idea for this method was based on the work that a student completed for his undergraduate research project in chemistry!

Because of the importance of liquid fuels in this country, fuels other than ethanol and biodiesel are also being researched. There is also considerable interest in the

use of plug-in electrical cars that would not require any liquid fuels. If the electrical energy used to charge the batteries in electric cars comes from renewable sources of electricity such as wind, solar, or geothermal, then the need for liquid fuels could be greatly decreased.


In 2007, the United States consumed a combined total of about 7.5 billion gallons of ethanol and biodiesel. By comparison, about 140 billion gallons of gasoline and 40 billion gallons of diesel fuel were consumed. Therefore, biofuels presently represent a small percentage of our total fuel consumption. Recently, Congress passed a bill requiring 36 billion gallons of biofuel to be produced yearly by 2022. Even if this goal is met, it is likely that we will still primarily rely on both fossil fuels and biofuel for the foreseeable future.

REFERENCES


- Biello, D. Grass Makes Better Ethanol than Corn Does. *Sci. Am.* [Online] **2008**, (Jan).
- Bomgardner, M. Seeking Better Biofuels. *Chem. Eng. News*, **2015**, 93 (Oct 26), 8
- Dale, B. E.; Pimentel, D. Point/Counterpoint: The Costs of Biofuels. *Chem. Eng. News* **2007**, 85 (Dec 17), 12.
- Fargione, J.; Hill, J.; Tilman, D.; Polasky, S.; Hawthorne, P. Land Clearing and Biofuel Carbon Debt. *Science* **2008**, 319 (Feb 29), 1235.
- Grunwald, M. The Clean Energy Scam. *Time* **2008**, 171 (Apr 7), 40.
- Heywood, J. B. Fueling Our Transportation Future. *Sci. Am.* **2006**, 295 (Sep), 60.
- Kammen, D. M. The Rise of Renewable Energy. *Sci. Am.* **2006**, 295 (Sep), 84.
- Kram, J. W. Minnesota Scientists Create New Biodiesel Manufacturing Process. *Biodiesel Magazine*, (Apr 7, 2008).
- Searchinger, T. Use of U.S. Croplands for Biofuels Increases Greenhouse Gases Through Emissions from Land-Use Change. *Science* **2008**, 319 (Feb 29), 1238–1248.

Biodiesel¹

In this experiment, you will prepare biodiesel from a vegetable oil in a base-catalyzed transesterification reaction:

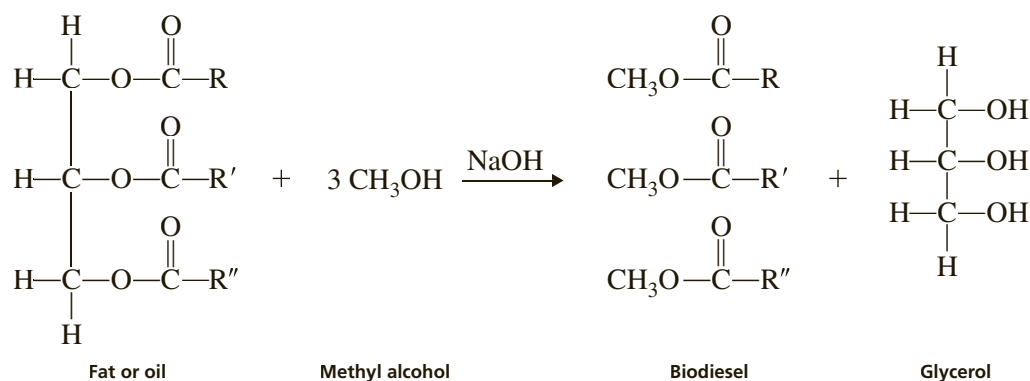

The first step in the mechanism for this synthesis is an acid-base reaction between sodium hydroxide and methyl alcohol:

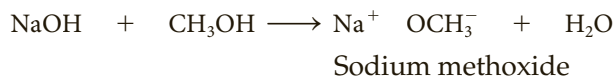


Methoxide ion is a strong nucleophile that now attacks the three carbonyl groups in the vegetable oil molecule. In the last step, glycerol and biodiesel are produced.

Because the R groups may have different numbers of carbons and they may be saturated (no double bonds) or may have one or two double bonds, biodiesel is a mixture of different molecules—all of which are methyl esters of fatty acids that made up the original vegetable oil. Most of the R groups have 10–18 carbons that are arranged in straight chains.

When the reaction is complete, the mixture is cooled and then centrifuged in order to separate the layers more completely. Since some unreacted methyl alcohol will be dissolved in the biodiesel layer, this layer is heated in an open container to remove all the methyl alcohol. The remaining liquid should be pure biodiesel.

¹ This experiment is based on a similar experiment developed by John Thompson, Lane Community College, Eugene, Oregon. It is posted on Greener Educational Materials (GEMs), an interactive database on green chemistry that is found on the University of Oregon green chemistry website (<http://greenchem.uoregon.edu/>).

When biodiesel is burned as a fuel, the following reaction occurs:

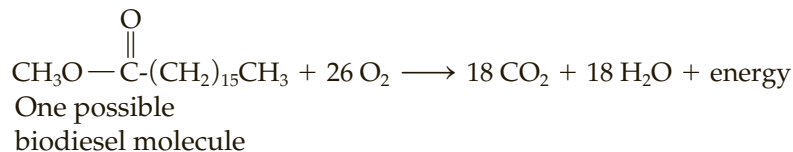

Burning biodiesel will produce a specific amount of energy, which can be measured using a bomb calorimeter. By combusting a specific weight of your biodiesel and measuring the temperature increase of the calorimeter, you can calculate the heat of combustion of biodiesel.

In Experiment 28A, coconut oil is converted into biodiesel, and other oils are converted into biodiesel in Experiment 28B. In Experiment 28C, the biodiesel is analyzed by infrared spectroscopy, NMR spectroscopy, and gas chromatography-mass spectrometry (GC-MS). The heat of combustion of biodiesel can also be determined in Experiment 28C.

REQUIRED READING

- New:* Technique 22 Gas Chromatography, Section 22.13
 Technique 25 Infrared Spectroscopy
 Technique 26 Nuclear Magnetic Resonance Spectroscopy
- Essays:* Biofuels
 Fats and Oils

SUGGESTED WASTE DISPOSAL

Discard the glycerol layer and leftover biodiesel into the container for the disposal of nonhalogenated organic waste.

NOTES TO THE INSTRUCTOR

We have found this experiment to be a good way to introduce infrared spectroscopy, NMR spectroscopy, and GC-MS. It is helpful to place the bottle containing the coconut oil into a beaker of warm water to keep the oil in the liquid state.

A Hewlett Packard model 5890 gas chromatograph-mass spectrometer may be used to obtain the GC-MS spectra of the biodiesel samples. Use a 30-meter Rtx-5 (Fused Silica) column (0.25 mm ID, 0.25 micron film). Set the inlet temperature at 250°C and the detector temperature at 280°C. The detection solvent delay is 3.3 minutes. The oven program has an initial temperature of 80°C (initial time 3 minutes) and a final temperature of 280°C (final time 3 minutes) with a rate of 20.0°C/minute. The sample is prepared by adding one drop of biodiesel to 4 mL of methanol and a 1.0 μL sample is injected. The GC-MS spectrum for biodiesel from coconut oil has several well-resolved peaks, whereas the spectra for biodiesel from other oils usually include overlapping peaks.

28A EXPERIMENT 28A*Biodiesel from Coconut Oil***PROCEDURE**

Prepare a warm-water bath in a 250-mL beaker. Use about 50 mL of water and heat the water to 55–60°C on a hot plate. (Do not let the temperature exceed 60° during the reaction period.) Weigh a 25-mL round-bottom flask. Add 10 mL of coconut oil to the flask and reweigh to get the weight of the oil. (Note: The coconut oil must be heated slightly in order to convert it to a liquid that can be measured in a graduated cylinder. It may also be advisable to warm the graduated cylinder.) Transfer 2.0 mL of sodium hydroxide dissolved in methyl alcohol solution to the flask.² (Note: Swirl the sodium hydroxide mixture before taking the 2-mL portion to make sure that the mixture is homogenous.) Place a magnetic stir bar in the round-bottom flask and attach the flask to a water condenser. (You do not need to run water through the water condenser.) Clamp the condenser so that the round-bottom flask is close to the bottom of the beaker. Turn on the magnetic stirrer to the highest level possible (this may not be the highest setting on the stirrer if the stir bar does not spin smoothly at high speeds). Stir for 30 minutes.

Transfer all of the liquid in the flask to a 15-mL plastic centrifuge tube with a cap and let it set for about 15 minutes. The mixture should separate into two layers: the larger top layer is biodiesel and the lower layer is mainly glycerol. To separate the layers more completely, place the tube in a centrifuge and spin for about 5 minutes (don't forget to counterbalance the centrifuge). If the layers have not separated completely after centrifugation, continue to centrifuge for another 5–10 minutes at a higher speed.

Using a Pasteur pipette, carefully remove the top layer of biodiesel and transfer this layer to a preweighed 50-mL beaker. You should leave behind a little of the biodiesel layer to make sure you don't contaminate it with the bottom layer.

Place the beaker on a hot plate and insert a thermometer into the bio-diesel, holding the thermometer in place with a clamp. Heat the biodiesel to about 70°C for 15–20 minutes to remove all the methyl alcohol. When the biodiesel has cooled to room temperature, weigh the beaker and liquid and calculate the weight of biodiesel produced. Record the appearance of the biodiesel.

To analyze your biodiesel, proceed to Experiment 28C.

² Note to instructor: Dry sodium hydroxide pellets overnight in an oven at 100°C. After grinding the dried sodium hydroxide with a mortar and pestle, add 0.875 g of this to an Erlenmeyer flask containing 50 mL of highly pure methanol. Place a magnetic stir bar in the flask and stir until all of the sodium hydroxide has dissolved. The mixture will be slightly cloudy.

28B EXPERIMENT 28B*Biodiesel from Other Oils*

Follow the procedure in Experiment 28A (Biodiesel from Coconut Oil), except use a different oil than coconut. Any of the oils listed at the bottom of Table 2 in the essay “Fats and Oils” that precedes Experiment 25 can be used. It will not be necessary to heat the oil when measuring out the 10 mL of oil, as all of these oils are liquids at room temperature.

To analyze your biodiesel, proceed to Experiment 28C.

28C EXPERIMENT 28C*Analysis of Biodiesel*

Spectroscopy. Obtain an infrared spectrum using salt plates (see Technique 25, Section 25.2). Determine the proton NMR spectrum using 3–4 drops of your biodiesel in 0.7 mL of deuterated chloroform. Since biodiesel consists of a mixture of different molecules, it is not helpful to perform an integration of the area under the peaks. Compare the NMR spectrum of biodiesel to the spectrum of vegetable oil shown here. Finally, analyze your sample using gas chromatography-mass spectrometry (GC-MS). Your instructor will provide you with instructions on how to do this.

Calorimetry (optional). Determine the heat of combustion (in kJoules/gram) of your biodiesel. Your instructor will provide instructions on how to use the bomb calorimeter and how to perform the calculations.


REPORT

Calculate the percent yield of biodiesel. This is difficult to do in the normal way based on moles because the vegetable oil and biodiesel molecules have variable composition. Therefore, you can base this calculation on the weight of oil used and the weight of biodiesel produced.

Analyze the infrared spectrum by identifying the principal absorption bands. Look for peaks in the spectrum that may indicate possible contamination from methanol, glycerol, or free fatty acids. Indicate any impurities found in your biodiesel bases on the infrared spectrum.

Analyze the NMR spectrum by comparing it to the NMR spectrum of vegetable oil with some of the signals labeled that is shown below. Look for evidence in the NMR spectrum for contamination by methanol, free fatty acids, or the original vegetable oil. Indicate any impurities found based on the NMR spectrum.

The library search contained in the software for the GC-MS instrument will give you a list of components detected in your sample, as well as the retention time and relative area (percentage) for each component. The results will also list possible substances that the computer has tried to match against the mass spectrum of each component. This list—often called a “hit list”—will include the name of each

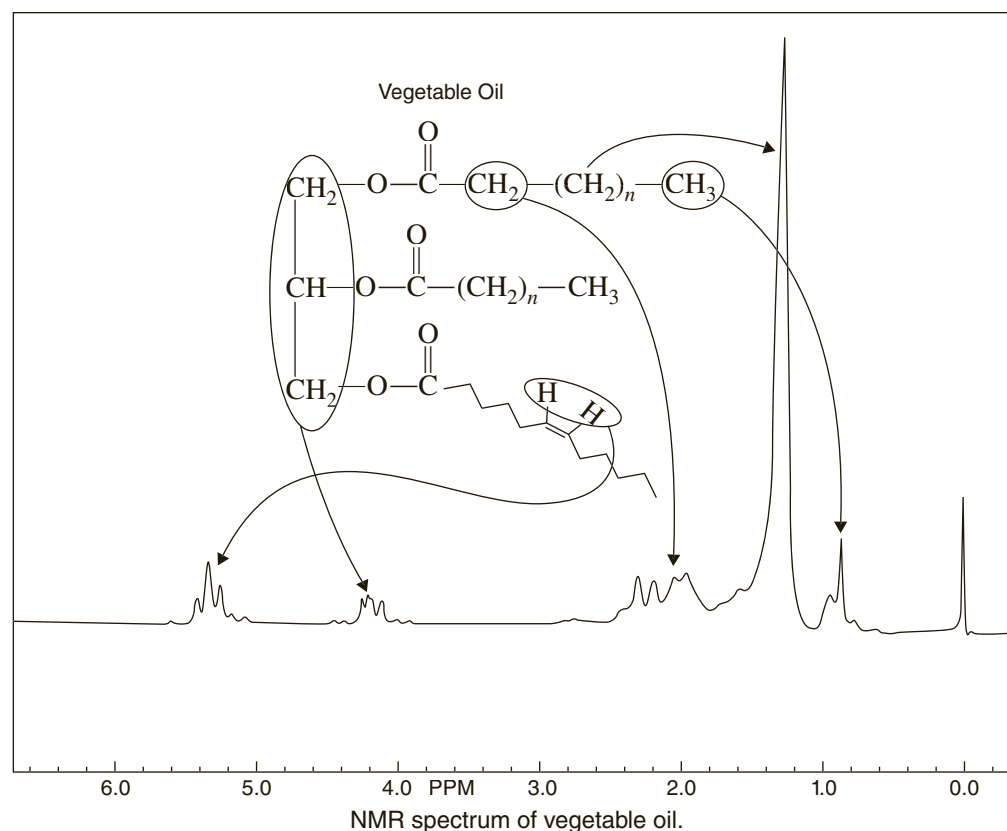


possible compound, its Chemical Abstracts Registry number (CAS number), and a “quality” (“confidence”) measure expressed as a percentage. The “quality” parameter estimates how closely the mass spectrum of the substance on the “hit list” fits the observed spectrum of that component in the gas chromatogram. The components that you identify from the GC-MS will be the methyl esters of the fatty acids that were initially part of the vegetable oil molecule. From the GC-MS data, you can determine the fatty acid composition (by percentages) in the original vegetable oil. Make a table of the main fatty acid components and the relative percentages. Compare this with the fatty acid composition given for this oil in Table 2 in the essay “Fats and Oils” that precedes Experiment 25. Is the fatty acid composition the same, and how do the relative percentages compare?

If you performed the experiment with the bomb calorimeter, list the data and calculate the heat of combustion for biodiesel in kJ/g. The heat of combustion for heptane, a component of gasoline, is 45 kJ/g. How do they compare?

QUESTIONS

- Write a complete reaction mechanism for this base-catalyzed transesterification reaction. Rather than starting with a complete oil molecule, give the mechanism for the reaction between the following ester and methanol in the presence of NaOH.

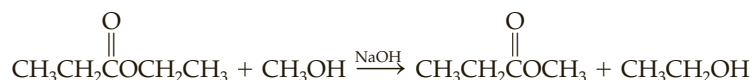

2. If you calculated the heat of combustion of biodiesel using a bomb calorimeter, answer the following question:

Compare the heat of combustion of biodiesel with heptane. Why does heptane have a larger heat of combustion? The heat of combustion of heptane is 45 kJ/g. (*Hint:* In answering this question, it may be helpful to compare the molecular formulas of biodiesel and heptane).

3. One argument for using biodiesel rather than gasoline is that the net amount of carbon dioxide released into the atmosphere from combusting biodiesel is sometimes claimed to be zero (or near zero). How can this argument be made, given that the combustion of biodiesel also releases carbon dioxide?
4. When the reaction for making biodiesel occurs, two layers are formed: biodiesel and glycerol. In which layer will most of each of the following substances be found? If a substance will be found to a large extent in both layers, you should indicate this.



Green Chemistry

The economic prosperity of the United States demands that it continue to have a robust chemical industry. In this age of environmental consciousness, however, we can no longer afford to allow the type of industry that has been characteristic of past practices to continue operating as it always has. There is a real need to develop an environmentally benign, or “green,” technology. Chemists must not only create new products but also design the chemical syntheses in a way that carefully considers their environmental ramifications.

Beginning with the first Earth Day celebration in 1970, scientists and the general public began to understand that the earth is a closed system in which the consumption of resources and indiscriminate disposal of waste materials are certain to bring about profound and long-lasting effects on the worldwide environment. Over the past decade, interest has begun to grow in an initiative known as *green chemistry*.

Green chemistry may be defined as the invention, design, and application of chemical products and processes to reduce or eliminate the use and generation of hazardous substances. Practitioners of green chemistry strive to protect the environment by cleaning up toxic waste sites and by inventing new chemical methods that do not pollute and that minimize the consumption of energy and natural resources. Guidelines for developing green chemistry technologies are summarized in the “Twelve Principles of Green Chemistry” shown in the table.

THE TWELVE PRINCIPLES OF GREEN CHEMISTRY

1. It is better to prevent waste than to treat or clean up waste after it is formed.
2. Synthetic methods should be designed to maximize the incorporation of all materials used in the process into the final product.
3. Wherever practicable, synthetic methodologies should be designed to use and generate substances that possess little or no toxicity to human health and the environment.
4. Chemical products should be designed to preserve efficacy of function while reducing toxicity.
5. The use of auxiliary substances (solvents, separation agents, etc.) should be made unnecessary whenever possible and innocuous when used.
6. Energy requirements should be recognized for their environmental and economic impacts and should be minimized. Synthetic methods should be conducted at ambient temperature and pressure.
7. A raw material or feedstock should be renewable rather than depleting whenever technically and economically practicable.
8. Unnecessary privatization (blocking group, protection/deprotection, temporary modification of physical/chemical processes) should be avoided whenever possible.
9. Catalytic reagents (as selective as possible) are superior to stoichiometric reagents.

10. Chemical products should be designed so that at the end of their function they do not persist in the environment and do break down into innocuous degradation products.
11. Analytical methodologies need to be further developed to allow for real-time, in-process monitoring and control before the formation of hazardous substances.
12. Substances and the form of a substance used in a chemical process should be chosen to minimize the potential for chemical accidents, including releases, explosions, and fires.

Source: Anastas, P. T.; Warner, J. C. *Green Chemistry: Theory and Practice*; Oxford University Press: New York, 1998. Reprinted by permission of the publisher.

The green chemistry program was begun shortly after the passage of the Pollution Prevention Act of 1990 and is the central focus of the Environmental Protection Agency's Design for the Environment Program. As a stimulus for research in the area of reducing the impact of chemical industry on the environment, the Presidential Green Chemistry Challenge Award was begun in 1995. The theme of the Green Chemistry Challenge is "Chemistry is not the problem; it's the solution." Since 1995, award winners have been responsible for the elimination of more than 460 million pounds of hazardous chemicals and have saved more than 440 million gallons of water and 26 million barrels of oil.

Winners of the Green Chemistry Challenge Award have developed foam fire retardants that do not use halons (compounds containing fluorine, chlorine, or bromine), cleaning agents that do not use tetrachloroethylene, methods that facilitate the recycling of polyethylene terephthalate soft-drink bottles, a method of synthesizing ibuprofen that minimizes the use of solvents and the generation of wastes, and a formulation that promotes the efficient release of ammonia from urea-based fertilizers. This latter contribution allows a more environmentally friendly means of applying fertilizers without the need for tilling or disturbing (and losing) precious topsoil.

Green syntheses of the future will require making choices about reactants, solvents, and reaction conditions that are designed to reduce resource consumption and waste production. We need to think about performing a synthesis in a way that will not consume excessive amounts of resources (and thus use less energy and be more economical), that will not produce excessive amounts of toxic or harmful by-products, and that will require milder reaction conditions.


The application of green-chemistry principles in an organic synthesis begins with the selection of the starting materials, called **feedstock**. Most organic compounds used as feedstock are derived from petroleum, a nonrenewable resource (see essay "Petroleum and Fossil Fuels" that precedes Experiment 27). A green approach is to replace these petrochemicals with chemicals derived from biological sources such as trees, corn, or soybeans. Not only is this approach more sustainable, but the refining of organic compounds from these plant-derived materials, sometimes called **biomass**, is also less polluting than the refining process for petrochemicals. Many pharmaceuticals, plastics, agricultural chemicals, and even transportation fuels can now be produced from chemicals derived from biomass. A good example of this is adipic acid, an organic chemical widely used in the production of nylon and lubricants. Adipic acid can be produced from benzene, a toxic petrochemical, or from glucose, which is found in plant sources.

Industrial processes are being designed that are based on the concept of **atom economy**. Atom economy means that close attention is paid to the design of chemical reactions so that all or most of the atoms that are starting materials in the process are converted into molecules of the desired product rather than into wasted by-products. Atom economy in the industrial world is the equivalent of ensuring that a chemical reaction proceeds with a high percentage yield in a classroom laboratory experiment.

The atom economy for a reaction can be calculated using the following equation:

$$\text{Percent atom economy} = \frac{\text{Molecular weight of desired product}}{\text{Molecular weights of all reactants}} \times 100\%$$

For example, consider the reaction for the synthesis of aspirin (Experiment 9, "Acetylsalicylic Acid"):

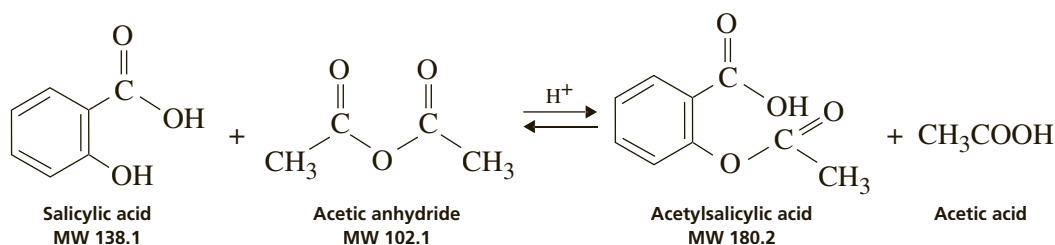

$$\text{Percent atom economy} = \frac{180.2}{138.1 + 102.1} \times 100 = 75.0$$

This calculation assumes the complete conversion of reactants into product and 100% recovery of the product, which is not possible. Furthermore, the calculation does not take into account that often an excess of one reactant is used to drive the reaction to completion. In this reaction, acetic anhydride is used in large excess to ensure the production of more acetylsalicylic acid. Nonetheless, the atom economy calculation is a good way to compare different possible pathways to a given product.


To illustrate the benefits of atom economy, consider the synthesis of ibuprofen, mentioned earlier, that won the Presidential Green Chemistry Challenge Award in 1997. In the former process, developed in the 1960s, only 40% of the reactant atoms were incorporated into the desired ibuprofen product; the remaining 60% of the reactant atoms found their way into unwanted by-products or wastes that required disposal. The new method requires fewer reaction steps and recovers 77% of the reactant atoms in the desired product. This "green" process eliminates millions of pounds of waste chemical by-products every year, and it reduces by millions of pounds the amount of reactants needed to prepare this widely used analgesic.

Another green chemistry approach is to select safer reagents that are used to carry out the synthesis of a given organic compound. In one example of this, milder or less toxic oxidizing reagents may be selected to carry out a conversion that is normally done in a less green way. For example, Oxone[®] can be used in some oxidation reactions instead of the highly toxic dichromate-sulfuric acid mixture. In some reactions, it is possible to use biological reagents, such as enzymes, to carry out a transformation. Another approach in green chemistry is to use a reagent that can promote the formation of a given product in less time and with greater yield. Finally, some reagents, especially catalysts, can be recovered at the end of the reaction period and recycled for use again in the same conversion.

Many solvents used in traditional organic syntheses are highly toxic. The green chemistry approach to the selection of solvents has resulted in several strategies. One method that has been developed is to use supercritical carbon dioxide as a solvent. Supercritical carbon dioxide is formed under conditions of high pressure in which the gas and liquid phases of carbon dioxide combine to a single-phase compressible fluid that becomes an environmentally benign solvent (temperature 31°C; pressure 7280 kPa, or 72 atm). Supercritical CO₂ has remarkable properties. It behaves as a material whose properties are intermediate between those of a solid and those of a liquid. The properties can be controlled by manipulating temperature and pressure. Supercritical CO₂ is environmentally benign because of its low toxicity and easy recyclability. Carbon dioxide is not added to the atmosphere; rather, it is removed from the atmosphere for use in chemical processes. It is used as a medium to carry out a large number of reactions that would otherwise have many negative environmental consequences. It is even possible to perform stereoselective synthesis in supercritical CO₂.

Some reactions can be carried out in ordinary water, the most green solvent possible. Recently, there has been much success in using near-critical water at higher temperatures where water behaves more like an organic solvent. Two of the award winners of the 2004 Green Chemistry Award, Charles Eckert and Charles Liotta, have advanced our understanding of supercritical CO₂ and near-critical water as solvents. One example of their work takes advantage of the dissociation of water that takes place under near-critical conditions, leading to a high concentration of hydronium and hydroxide ions. These ions can serve as self-neutralizing catalysts, and they can replace catalysts that must normally be added to the reaction mixture. Eckert and Liotta were able to run Friedel-Crafts reactions (Experiment 60, “Friedel-Crafts Acylation”) in near-critical water without the need for the acid catalyst AlCl₃, which is normally used in large amounts in these reactions.

Research has also focused on **ionic liquids**, salts that are liquid at room temperature and do not evaporate. Ionic liquids are excellent solvents for many materials, and they can be recycled. An example of an ionic liquid is

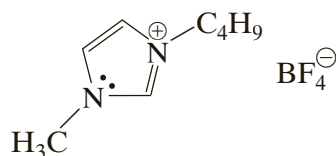

Even though many of the ionic liquids are expensive, their high initial cost is mitigated because, through recycling, they are not consumed or discarded. In addition, product recovery is often easier than with traditional solvents. In the past five years, many new ionic liquids have been developed with a broad range of properties. By selecting the appropriate ionic liquid, it is now possible to carry out many types of organic reactions in these solvents. In some reactions, a well-designed ionic solvent can lead to better yields under milder conditions than is possible with traditional solvents. Recently, researchers have developed ionic liquids made from artificial “sweeteners” that are nontoxic and extend even further the concept of green chemistry.

It is possible in some organic syntheses to completely eliminate the need for any solvent! Some reactions that are traditionally carried out in solvents can be carried out either in the solid or gas phases without the presence of any solvent.

Another approach to making organic chemistry greener involves the way in which a reaction is carried out, rather than in the selection of starting material, reagents, or solvents. Microwave technology (see Technique 7.12, Microwave-Assisted

Organic Chemistry) can be used in some reactions to provide the heat energy required to make the transformation go to completion. With microwave technology, reactions can take place with less-toxic reagents, in a shorter time, and with fewer side reactions—all goals of green chemistry. Microwave technology has also been used to create supercritical water that behaves more like an organic solvent and could replace more toxic solvents in carrying out organic reactions.

Another green approach involving technology is the use of solid-phase extraction (SPE) columns (see Technique 12.14, Solid Phase Extraction). Using SPE columns, extractions such as removing caffeine from tea can be carried out more quickly and with less-toxic solvents. In other applications, SPE columns can be used to carry out the synthesis of organic compounds more efficiently with less use of toxic reagents.

Industry has discovered that environmental stewardship makes good economic sense, and there is a renewed interest in cleaning up manufacturing processes and products. In spite of the continuing adversarial nature of relations between industry and environmentalists, companies are discovering that preventing pollution in the first place, using less energy, and developing atom-economic methods makes as much sense as spending less money on raw materials or capturing a greater share of the market for their product. Although U.S. chemical industries are by no means near their stated goal of reducing the emission of toxic substances to zero or near-zero levels, significant progress is being made.

The teaching of the principles of green chemistry is beginning to find its way into the classroom. In this textbook, we have attempted to improve the green qualities of some of the experiments and have added several green experiments. The following table lists the experiments in this textbook that have a significant green component, along with the primary aspect of the experiment that makes it green.

<i>Experiment</i>	<i>Green Aspect</i>
Exp. 27, "Gas Chromatographic Analysis of Gasolines"	Discussion of pollution-controlling additives
Exp. 28, "Biodiesel"	Transportation fuel using recycled materials
Exp. 29, "Chiral Reduction of Ethyl Acetoacetate"	Biological reagent, baker's yeast
Exp. 30, Nitration of Aromatic Compounds Using a Recyclable Catalyst"	Use of a recyclable catalyst to increase reaction efficiency
Exp. 31, "Reduction of Ketones Using Carrot Extract"	Biological reagent
Exp. 33, "An Oxidation–Reduction Scheme: Borneol, Camphor, Isoborneol"	Less-toxic oxidizing agents
Exp. 36, "Aqueous-Based Organozinc Reactions"	Water used as the solvent
Exp. 37, "Sonogashira Coupling of Iodoaromatic Compounds with Alkynes"	Use of a recyclable catalyst to increase reaction efficiency
Exp. 38, "Grubb's-Catalyzed Metathesis of Eugenol with <i>cis</i> -1,4-Butenediol"	Use of a recyclable catalyst to increase reaction efficiency
Exp. 40, "A Green Enantioselective Aldol Condensation Reaction"	Use of less-toxic reagents
Exp. 42, "Preparation of Triarylpyridines"	Solvent-less reaction

<i>Experiment</i>	<i>Green Aspect</i>
Exp. 43, "1,4-Diphenyl-1,3-butadiene"	Solvent-less reaction
Exp. 51, "Synthesis of a New Polymer Using Grubb's Catalyst"	Use of a recyclable catalyst to increase reaction efficiency
Exp. 54, "Diels-Alder Reaction with Anthracene-9-methanol"	Water used as the solvent
Exp. 66, "Green Epoxidation of Chalcones"	Use of a less-toxic reagents
Exp. 67, "Cyclopropanation of Chalcones"	Use of a less-toxic reagents

In addition, Experiment 57 (Identification of Unknowns) offers a "green" alternative procedure. This procedure avoids the use of toxic chemicals for classification tests and substitutes the use of spectroscopy, which does not require any chemical reagents (except a small amount of organic solvent).

Certainly, enormous challenges remain. Generations of new scientists must be taught that it is important to consider the environmental impact of any new methods that are introduced. Industry and business leaders must learn to appreciate that adopting an atom-economic approach to the development of chemical processes makes good long-term economic sense and is a responsible means of conducting business. Political leaders must also develop an understanding of what the benefits of a green technology can be and why it is responsible to encourage such initiatives.

REFERENCES

- Amato, I. Green Chemistry Proves It Pays Companies to Find New Ways to Show That Preventing Pollution Makes More Sense Than Cleaning Up Afterward. *Fortune* [Online], Jul 24, 2000. www.fortune.com/fortune/articles/0.15114.368198.00.html
- Freemantle, M. Ionic Liquids in Organic Synthesis. *Chem. Eng. News* **2004**, 82 (Nov 8), 44.
- Jacoby, M. Making Olefins from Soybeans. *Chem. Eng. News* **2005**, 83 (Jan 3), 10.
- Mark, V. Riding the Microwave. *Chem. Eng. News* **2004**, 82 (Dec 13), 14.
- Matlack, A. Some Recent Trends and Problems in Green Chemistry. *Green Chem.* **2003** (Feb), G7-G11.
- Mullin, R. Sustainable Specialties. *Chem. Eng. News* **2004**, 82 (Nov 8), 29.
- Oakes, R. S.; Clifford, A. A.; Bartle, K. D.; Pett, M. T.; Rayner, C. M. Sulfur Oxidation in Supercritical Carbon Dioxide: Dramatic Pressure Dependent Enhancement of Diastereo-selectivity for Sulfoxidation of Cysteine Derivatives. *Chem. Comm.* **1999**, 247-248.
- Ritter, S. K. Green Innovations. *Chem. Eng. News* **2004**, 82 (Jul 12), 25.

Chiral Reduction of Ethyl Acetoacetate; Optical Purity Determination Using a Chiral Shift Reagent

Green chemistry

Stereochemistry

Reduction with yeast

Use of a separatory funnel

Chiral gas chromatography


Polarimetry

Optical purity (enantiomeric excess) determination


Nuclear magnetic resonance (optional)

Chiral chemical shift reagents (optional)

The experiment described in Experiment 29A uses common baker's yeast as a chiral reducing medium to transform an achiral starting material, ethyl acetoacetate, into a chiral product. When a single stereoisomer is formed in a chemical reaction from an achiral starting material, the process is said to be **enantiospecific**. In other words, one stereoisomer (enantiomer) is formed in preference to its mirror image. In the present experiment, the ethyl (*S*)-3-hydroxybutanoate stereoisomer is formed preferentially. In actual fact, however, some of the (*R*)-enantiomer is also formed in the reaction. The reaction, therefore, is described as an **enantioselective** process because the reaction does not produce one stereoisomer exclusively. Chiral gas chromatography and polarimetry will be employed to determine the percentages of each of the enantiomers. Generally, the chiral reduction produces less than 8% of the ethyl (*R*)-3-hydroxybutanoate.

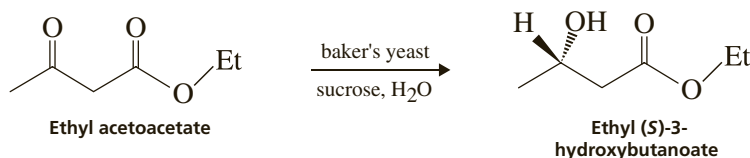

In contrast, when ethyl acetoacetate is reduced with sodium borohydride in methanol, the reaction yields a 50–50 mixture of the (*R*) and (*S*)-stereoisomers. A racemic mixture is formed because the reaction is not being conducted in a chiral medium.

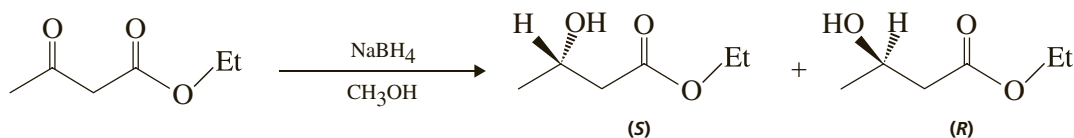

We are grateful to Dr. Snorri Sigurdsson and James Patterson, University of Washington, Seattle, for suggested improvements.

In Experiment 29B (optional), you may use nuclear magnetic resonance spectroscopy to determine the relative amounts of (*R*) and (*S*) enantiomers produced in the chiral reduction of ethyl acetoacetate. This part of the experiment requires the use of a chiral shift reagent.

29A EXPERIMENT 29A

Chiral Reduction of Ethyl Acetoacetate

REQUIRED READING

<i>Review:</i>	Technique 8	Filtration, Sections 8.3 and 8.4
	Technique 12	Extractions, Separations, and Drying Agents, Sections 12.4 and 12.10
	Techniques 22 and 25	
	Techniques 26 and 27 (optional)	
<i>New:</i>	Technique 23	Polarimetry
	Essay	Green Chemistry

SPECIAL INSTRUCTIONS

Day 1 of the experiment involves setting up the reaction. Another experiment can be conducted concurrently with this experiment. Part of this first laboratory period is used to mix the yeast, sucrose, and ethyl acetoacetate in a 500-mL Erlenmeyer flask. The mixture is stirred during part of that first period. The mixture is then covered and stored until the next period. The reduction requires at least 2 days.

Day 2 of the experiment is used to isolate the chiral ethyl 3-hydroxybutanoate. After this has been isolated, each student's product is analyzed by chiral gas chromatography and polarimetry to determine the percentages of each of the enantiomers. As an optional experiment (Experiment 29B), the products can also be analyzed by NMR using a chiral shift reagent to determine the percentages of each of the enantiomers present in the ethyl 3-hydroxybutanoate produced in the chiral reduction.

SUGGESTED WASTE DISPOSAL

The Celite, residual yeast, and cheesecloth from the reduction can be disposed of in the trash. The aqueous solutions and emulsion left from the extraction with methylene chloride should be placed in the aqueous waste container. Methylene chloride waste should be poured into the waste container designated for halogenated waste.

NOTES TO THE INSTRUCTOR

It is strongly advised that rotary evaporators be made available for this experiment. Approximately 90 mL of methylene chloride is used for each student. The experiment will be more "green" if the solvent can be recovered. The instructor

will need to make available to each student a large Büchner funnel (10 cm), 500-mL filter flask, 500-mL Erlenmeyer flask, 1.5- or 2-inch magnetic stir bar, and a 500-mL separatory funnel. It is advised that packaged dry yeast be used. We suggest Fleischmann's Rapid Rise (baker's) Yeast, which contains 7 g of yeast per package. Purchase packages of 100% cotton cheesecloth that consists of three layers (do not separate the layers from each other). Cut the three-ply cheesecloth into 4×8 inch strips to be folded into 4×4 inch sections for the Büchner funnel. In some cases, the yeast does not grow substantially during the first half hour. It is best to discard the mixture and start the reaction again if it appears that the yeast is not growing. In most cases, the temperature may not have been controlled carefully. It is recommended that the flasks containing the reaction mixture be stored in an area where the temperature is maintained at about 25°C , if possible. The optimal reduction period is four days. A small amount of unreduced ethyl acetoacetate remains after a 2-day reduction (less than 1%), and a 4-day reduction yields no remaining ethyl acetoacetate. The expected yield of chiral hydroxyester should be around 65%, consisting of 92–94% ethyl (*S*)-3-hydroxybutanoate.

PROCEDURE

Yeast Reduction

To a 500-mL Erlenmeyer flask, add 150 mL of de-ionized (DI) water and a 1.5- or 2-inch stir bar. Warm the water to about $35\text{--}40^{\circ}\text{C}$ using a hot plate set on low. Add 7 g of sucrose and 7 g of Fleischmann's Rapid Rise (dry baker's) Yeast to the flask. Swirl the contents of the flask in order to distribute the yeast into the aqueous solution; otherwise it will remain at the top of the solution. Stir the mixture for 15 minutes while maintaining the temperature at 35°C . During this time, the yeast will become activated and will grow substantially. Add 3.0 g of ethyl acetoacetate and 8 mL of hexane to the yeast mixture. Stir the mixture with a magnetic stirrer for 1.5 hours. Because the mixture may become thick, check periodically to see whether the mixture is being stirred. The reaction is somewhat exothermic, so you may not need to heat the mixture. Nevertheless, you should monitor the temperature to make sure that it remains near 30°C . Adjust the temperature to 30°C if the temperature falls below this value.

Label the Erlenmeyer flask with your name and ask your instructor to store the flask. Cover the top of the flask loosely with aluminum foil so that carbon dioxide can be expelled during the reduction. The mixture will stand, without stirring, until the next laboratory period (2–4 days). At some point during the laboratory period, obtain the infrared spectrum of ethyl acetoacetate for the purpose of comparison to the reduced product.

CAUTION



Do not breathe the Celite powder.

Isolation of the Alcohol Product

Obtain a 500-mL separatory funnel, a large Büchner (10 cm) funnel, and a 500-mL filter flask from your instructor. To the yeast solution, add 5 g of Celite and stir the mixture magnetically for 1 minute (see Technique 8, Section 8.4). Allow the solid to settle as much as possible (at least 5 minutes). Set up a vacuum-filtration apparatus using the large Büchner funnel (see Technique 8, Section 8.3). Wet one

piece of filter paper with water and place it into the funnel. Obtain a 4 × 8 inch strip of cheesecloth and fold it over to make a 4 × 4 inch square. Wet it with water and place it on top of the filter paper so that it completely covers the filter paper and is partly up the side of the Büchner funnel. You are now ready to filter your solution. Turn on the vacuum source (aspirator or the house vacuum). Decant the clear supernatant liquid slowly into the Büchner funnel. If you do this slowly, you may avoid plugging the filter paper with small particles. Once the supernatant liquid has been poured into the funnel, add the Celite slurry to the Büchner funnel. Rinse the flask with 20 mL of water and pour the remaining Celite–yeast mixture into the Büchner funnel. Discard the Celite, yeast, and cheesecloth waste into the trash. The Celite helps to trap the very tiny yeast particles. Some of the yeast and Celite will pass through the filter into the filter flask. This is unavoidable.

Add 20 g of sodium chloride to the filtrate in the filter flask and swirl the flask gently until the sodium chloride dissolves. If an emulsion forms, you may be swirling the flask too vigorously. Pour the filtrate into a 500-mL separatory funnel. Add 30 mL of methylene chloride to the funnel and stopper the funnel (see Technique 12, Section 12.4). In order to avoid a difficult emulsion, do not shake the separatory funnel; instead, slowly invert the funnel and bring it back to the upright position. Repeat this motion over a period of 5 minutes. Vent the funnel occasionally to relieve pressure. Drain the lower methylene chloride layer from the separatory funnel into a 250-mL Erlenmeyer flask, leaving behind a small amount of emulsion and the aqueous layer in the separatory funnel. Add another 30-mL portion of methylene chloride to the separatory funnel and repeat the extraction procedure. Drain the lower methylene chloride layer into the same Erlenmeyer flask holding the first methylene chloride extract. Repeat the extraction process a third time with a 30-mL portion of methylene chloride. Discard the emulsion and aqueous layer remaining in the separatory funnel into a suitable aqueous waste container.

Dry the three combined methylene chloride extracts over about 1 g of anhydrous granular sodium sulfate for at least 5 minutes. Occasionally, swirl the contents of the flask to help dry the solution. Decant the liquid into a 250-mL beaker and evaporate the solvent using an air or nitrogen stream until the volume of liquid remains constant (approximately 1–2 mL). (Alternatively, a rotary evaporator or distillation may be used to remove the methylene chloride from the product.)¹ Often the remaining liquid contains some water. To remove the water, add 10 mL of methylene chloride to dissolve the product and add 0.5 g of anhydrous granular sodium sulfate to the solution. Decant the methylene chloride solution away from the drying agent into a preweighed 50-mL beaker. Evaporate the solvent using an air or nitrogen stream until the volume of liquid remains constant. The liquid contains the ethyl (*S*)-3-hydroxybutanoate that has been produced by chiral reduction of ethyl acetoacetate. A small amount of ethyl acetoacetate may remain unreduced in the sample. Reweigh the beaker to determine the weight of the product. Calculate the percentage yield of product.

Infrared Spectroscopy

Determine the spectrum of your isolated product. The infrared spectrum provides the best direct evidence for the reduction of ethyl acetoacetate. Look for presence of a hydroxyl group (about 3440 cm⁻¹) that was produced in the reduction of the

¹ Pour the dry methylene chloride extracts into a round-bottom flask and remove the solvent with a rotary evaporator or by distillation. After removing the solvent, add 10 mL of fresh methylene chloride and 0.5 g of anhydrous granular sodium sulfate to the round-bottom flask. Decant the solution away from the drying agent into a preweighed beaker as indicated in the procedure.

carbonyl group. Compare the spectrum of the product, ethyl 3-hydroxybutanoate, to the starting material, ethyl acetoacetate. What differences do you notice in the two spectra? Label the two spectra with peak assignments and include them with your laboratory report.

Chiral Gas Chromatography

Chiral gas chromatography will provide a direct measure of amounts of each stereoisomer present in your chiral ethyl 3-hydroxybutanoate sample. A Varian CP-3800 equipped with an Alltech Cyclosil B capillary column (30 m, 0.25-mm ID, 0.25 μm) provides an excellent separation of (*R*) and (*S*)-enantiomers. Set the FID detector at 270°C and the injector temperature at 250°C, with a 50:1 split ratio. Set the column oven temperature at 90°C and hold at that temperature for 20 minutes. The helium flow rate is 1 mL/min. The compounds elute in the following order: ethyl (*S*)-3-hydroxybutanoate (14.3 min) and the (*R*)-enantiomer (15.0 min). Any remaining ethyl acetoacetate present in the sample will produce a peak with a retention time of 14.1 minutes. Your observed retention times may vary from those given here, but the order of elution will be the same. Calculate the percentages of each of the enantiomers from the chiral gas chromatography results. Usually, about 92–94% of the (*S*)-enantiomer is obtained from the reduction.

Polarimetry

Fill a 0.5-dm polarimeter cell with your chiral hydroxyester (about 2 mL required). You may need to combine your product with material obtained by one other student in order to have enough material to fill the cell. Determine the observed optical rotation for the chiral material. Your instructor will show you how to use the polarimeter. Calculate the specific rotation for your sample using the equation provided in Technique 23. The concentration value, *c*, in the equation is 1.02 g/mL. Using the published value for the *specific* rotation of ethyl (*S*)-(+)-3-hydroxybutanoate of $[\alpha_D^{25}] = +43.50$, calculate the optical purity (enantiomeric excess) for your sample (see Technique 23, Section 23.5). Report the observed rotation, the calculated specific rotation, the optical purity (enantiomeric excess), and the percentages of each of the enantiomers to the instructor. How do the percentages of each of the enantiomers calculated from the polarimeter measurement compare to the values obtained from chiral gas chromatography?²

Proton and Carbon NMR Spectroscopy (Optional)

At the option of the instructor, you may obtain the proton spectrum (shown in Figures 1 and 2 and interpreted in Experiment 29B) and the carbon NMR spectra of the product. The carbon NMR spectrum shows peaks at 14.3, 22.6, 43.1, 60.7, 64.3, and 172.7 ppm.

² The percentages calculated from polarimetry may vary considerably from those obtained by chiral gas chromatography. Often the samples contain some solvent and other impurities that reduce the observed optical rotation value. The solvent and impurities do not influence the more accurate percentages obtained directly by chiral gas chromatography.

29B EXPERIMENT 29B (OPTIONAL)

NMR Determination of the Optical Purity of Ethyl (S)-3-Hydroxybutanoate

In Experiment 29A, the yeast reduction of ethyl acetoacetate forms a product that is predominantly the (S)-enantiomer of ethyl 3-hydroxybutanoate. In this part of the experiment, we will use NMR to determine the percentages of each of the enantiomers in the product. The 300 MHz proton NMR spectrum of racemic ethyl 3-hydroxybutanoate is shown in Figure 1. The expansions of the individual patterns from Figure 1 are shown in Figure 2. The methyl protons (H_a) appear as a doublet at 1.23 ppm, and the methyl protons (H_b) appear as a triplet at 1.28 ppm. The methylene protons (H_c and H_d) are diastereotopic (nonequivalent) and appear at 2.40 and 2.49 ppm (each a doublet of doublets). The hydroxyl group appears at about 3.1 ppm. The quartet at 4.17 ppm results from the methylene protons (H_e) split by the protons (H_b). The methine proton (H_f) is buried under the quartet at about 4.2 ppm.

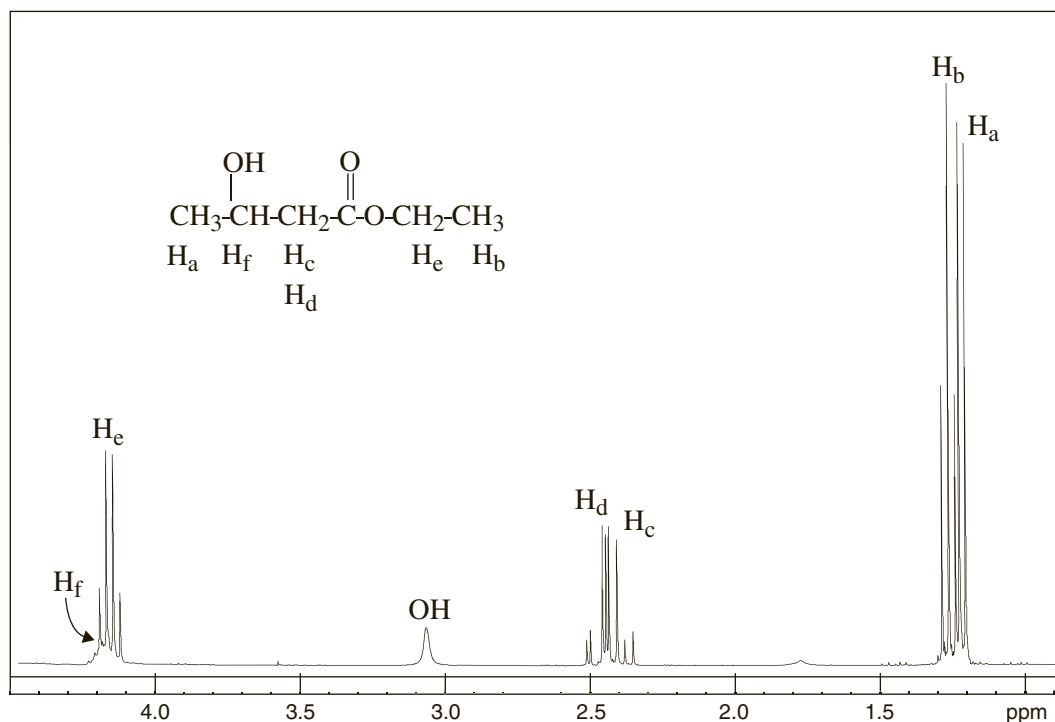
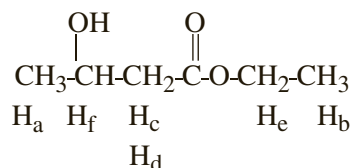


Figure 1
NMR spectrum (300 MHz) of racemic ethyl 3-hydroxybutanoate with no chiral shift reagent present.


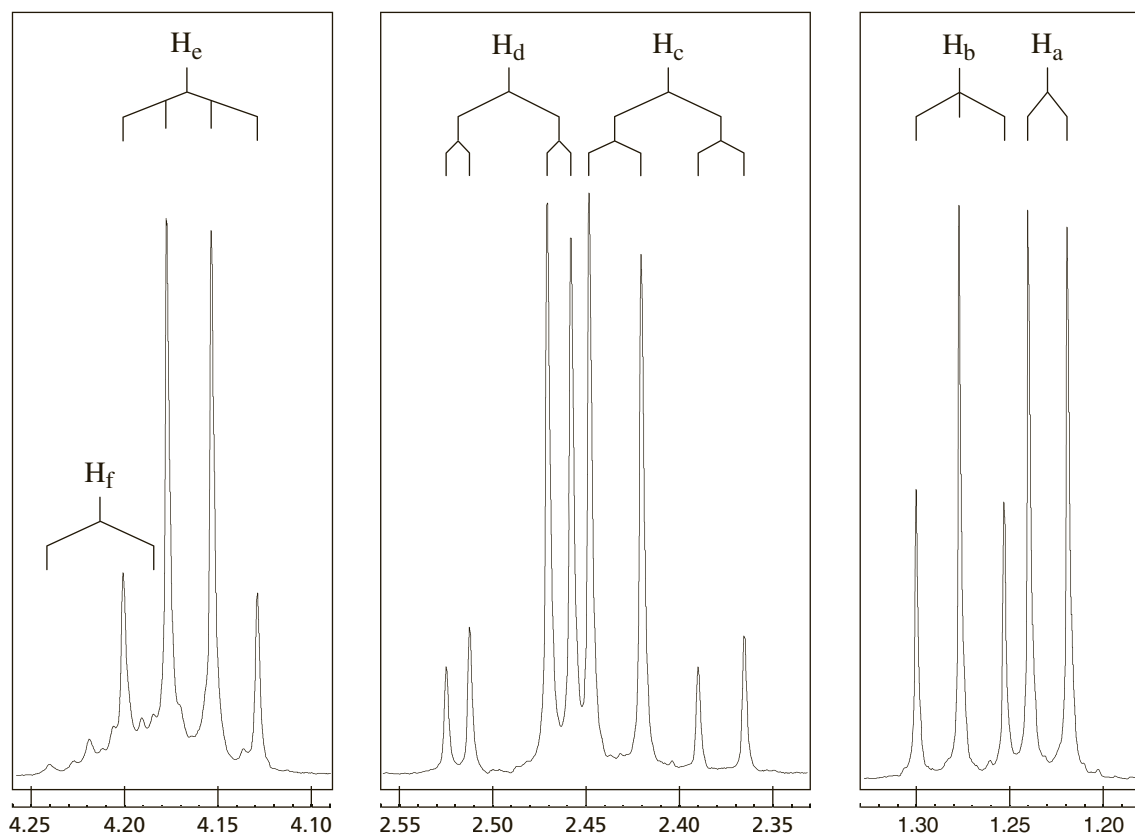


Figure 2
Expansions of the NMR spectrum of racemic ethyl 3-hydroxybutanoate.

Although the normal proton NMR spectrum for the racemic ethyl 3-hydroxybutanoate is not expected to be any different from the proton NMR spectra of each of the enantiomers in an **achiral environment**, the introduction of a chiral shift reagent creates a **chiral environment**. This chiral environment allows the two enantiomers to be distinguished from each other. A general discussion of nonchiral chemical shift reagents is found in Technique 26, Section 26.15. These reagents spread out the resonances of the compound with which they are used, increasing by the largest amount the chemical shifts of the protons that are nearest the center of the metal complex. Because the spectra of both enantiomers are identical under these conditions, the usual chemical shift reagent would not help our analysis. However, if we use a chemical shift reagent that is itself chiral, we can distinguish the two enantiomers by their NMR spectra. The two enantiomers, which are each chiral, will interact differently with the chiral shift reagent. The complexes formed from the (*R*) and (*S*)-enantiomers and with the chiral shift reagent will be diastereomers. Diastereomers usually have different physical properties, and the NMR spectra are no exception. The two complexes will be formed with slightly differing geometries. Although the effect may be small, it is large enough to see differences in the NMR spectra of the two enantiomers.

The chiral shift reagent used in this experiment is *tris*-[3-(heptafluoropropylhydroxymethylene)-(+)-camphorato]europium (III), or $\text{Eu}(\text{hfc})_3$. In this

complex, the europium is in a chiral environment because it is complexed to camphor, which is a chiral molecule. $\text{Eu}(\text{hfc})_3$ has the structure shown below.

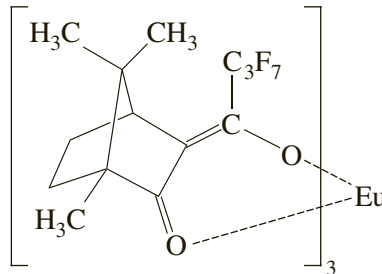

REQUIRED READING

New: Technique 26 Nuclear Magnetic Resonance Spectroscopy, Section 26.15

SPECIAL INSTRUCTIONS

This experiment requires the use of a high-field NMR spectrometer in order to obtain sufficient separation of peaks for the two enantiomers. The chiral shift reagent does cause some peak broadening, so care should be taken not to add too much of this reagent to the chiral ethyl 3-hydroxybutanoate sample. A 0.035-g sample of the chiral material and 8–11 mg of chiral shift reagent should be sufficient to give good results.


SUGGESTED WASTE DISPOSAL

Discard the remaining solution from your NMR tube into the container designated for the disposal of halogenated organic waste.

PROCEDURE

Using a Pasteur pipette to aid the transfer, weigh 0.035 g of chiral ethyl 3-hydroxybutanoate from Experiment 29A directly into an NMR tube. Weigh 8–11 mg of *tris* [3-(heptafluoropropylhydroxymethylene)-(+)-camphorato]europium(III) chiral shift reagent on a piece of weighing paper and add the chiral shift reagent to the chiral hydroxyester in the NMR tube. Take care to avoid chipping the fragile NMR tube while adding the shift reagent with a microspatula. Add CDCl_3 solvent to the NMR tube until the level reaches 50 mm. Cap the tube and invert it to mix the sample. Allow the NMR sample to stand for a minimum of about 5–8 minutes before determining the NMR spectrum. Record in your notebook the exact weights of sample and chiral shift reagent that you used.

Determine the NMR spectrum of the sample. The peaks of interest are the methyl protons, H_a (doublet) and H_b (triplet). Notice in Figure 3 that the doublet and triplet peaks for the two methyl groups in the *racemic* ethyl 3-hydroxybutanoate

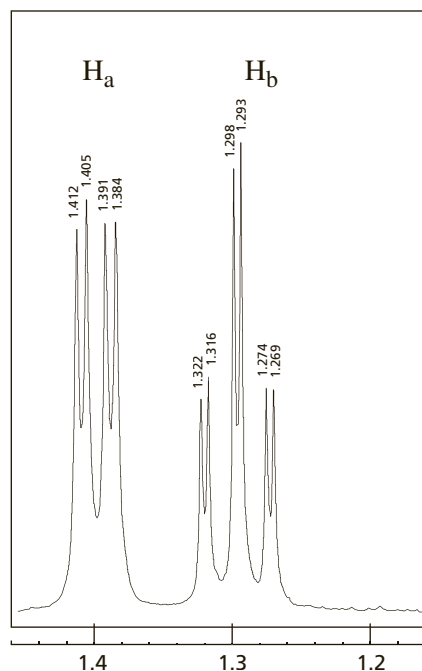
**Figure 3**

NMR spectrum (300 MHz) of racemic ethyl 3-hydroxybutanoate, with chiral shift reagent added.

Note: H_a for the (*S*)-enantiomer = 1.412, 1.391; H_b for the (*S*)-enantiomer = 1.322, 1.298, 1.274; H_a for the (*R*)-enantiomer = 1.405, 1.384; H_b for the (*R*)-enantiomer = 1.316, 1.293, 1.269.

are doubled. The downfield doublet (1.412 and 1.391 ppm) and triplet (1.322, 1.298, and 1.274 ppm) peaks are assigned to the (*S*)-enantiomer. The upfield doublet (1.405 and 1.384 ppm) and triplet (1.316, 1.293, and 1.269 ppm) peaks are assigned to the (*R*)-enantiomer. Your expansion of this area of the NMR spectrum should also show a doubling of the peaks as in Figure 3, but the upfield doublet for the (*R*)-enantiomer will be smaller. The same will be true for the (*R*)-enantiomer in the triplet pattern. By integration, determine the percentages of the (*S*)- and (*R*)-enantiomers in the chiral ethyl 3-hydroxybutanoate from Experiment 29A. Although the positions of the peaks may vary somewhat from those shown in Figure 3, you should still find that the doublet and triplet for the (*S*)-enantiomer will always be downfield relative to the (*R*)-enantiomer.

The assignments for the (*S*)- and (*R*)-enantiomers shown in Figure 3 were determined by obtaining the NMR spectrum of pure samples of each enantiomer in the presence of the chiral shift reagent (Figures 4 and 5). You may have noticed that the doublet has moved further downfield relative to the triplet (compare Figures 2 and 3). The reason for this is that the complexation of the chiral shift reagent occurs at the hydroxyl group. Because the methyl group (H_a) is closer to the europium

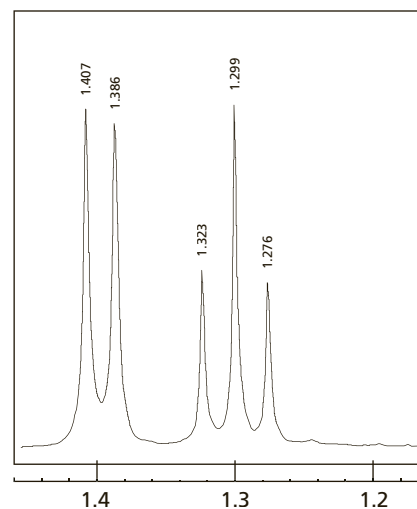**Figure 4**

NMR spectrum (300 MHz) of ethyl (*S*)-3-hydroxybutanoate, with chiral shift reagent added.

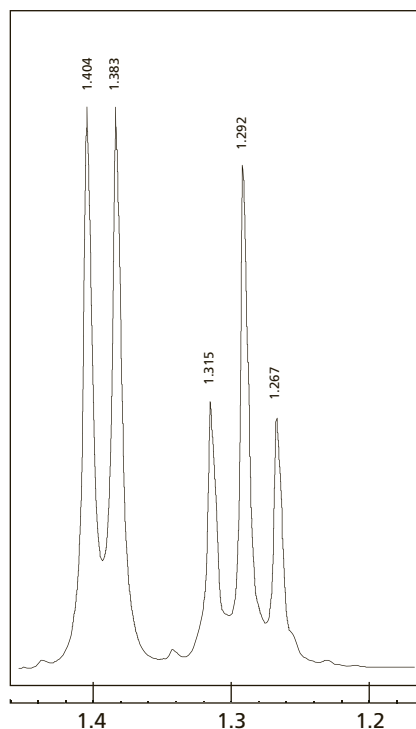


Figure 5
NMR spectrum (300 MHz) of ethyl (R)-3-hydroxybutanoate, with chiral shift reagent added.

atom, it is expected that that group will be shifted further downfield relative to the other methyl group (H_b).

REFERENCES

- Cui, J.-N.; Ema, T.; Sakai, T.; Utaka, M. Control of Enantioselectivity in the Baker's Yeast Asymmetric Reduction of Chlorodiketones to Chloro (S)-Hydroxyketones. *Tetrahedron: Asymmetry* **1998**, *9*, 2681–2691.
- Naoshima, Y.; Maeda, J.; Munakata, Y. Control of the Enantioselectivity of Bioreduction with Immobilized Bakers' Yeast in a Hexane Solvent System. *J. Chem. Soc. Perkin Trans. 1* **1992**, 659–660.
- Seebach, D.; Sutter, M. A.; Weber, R. H.; Züger, M. F. Yeast Reduction of Ethyl Acetoacetate: (S)-(+)-Ethyl 3-Hydroxybutanoate. *Org. Synth.* **1984**, *63*, 1–9.

QUESTIONS

1. Would you expect to see a difference in retention times for the ethyl (S)-3-hydroxybutanoate and the (R)-enantiomer on gas chromatography columns described in Technique 22?
2. What is the biological reducing agent that gives rise to the formation of chiral ethyl 3-hydroxybutanoate? You may need to look in a reference book to find an answer to this question.
3. Explain the NMR patterns for protons H_c and H_d shown in Figure 2. (*Hints:* These protons are nonequivalent because of their location adjacent to a stereocenter. The 2J coupling constants for protons attached to an sp^3 carbon are very large—in this case, 16.5 Hz. The 3J coupling constants are not equal. Draw a sawhorse projection for the molecule. Can you see why the 3J coupling constants might be different?)

Nitration of Aromatic Compounds Using a Recyclable Catalyst

Green chemistry

Nitration

Atom-economic reaction

Recyclable catalyst


Rotary evaporator (optional)

Mass spectrometry

Gas chromatography

Chemists in academia and industry are attempting to make chemical reactions more environmentally friendly (see the essay "Green Chemistry"). One way to accomplish this is to use exact (stoichiometric) amounts of starting reagents so that no excess material need be thrown away, thus contributing to a higher atom economy. Another aspect of Green Chemistry is that chemists should use catalysts. These materials have the advantage of allowing reactions to occur under milder conditions, and catalysts can also be reused. Thus, Green Chemistry helps keep the environment clean while producing useful products.

In the present experiment, we employ a Lewis acid, ytterbium (III) trifluoromethanesulfonate, as a catalyst for the nitration of a series of aromatic substrates with nitric acid. This catalyst will be recycled (recovered) and reused.



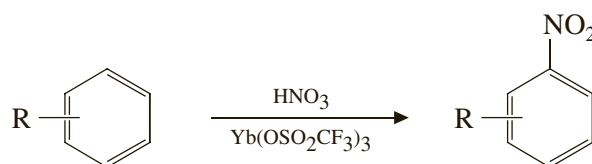


The solvent used in this reaction, 1,2-dichloroethane, is not environmentally friendly, but the solvent can be recovered using a rotary evaporator.

A proposed mechanism for this reaction involves the following three steps to generate the nitronium ion.¹ The trifluoromethanesulfonate (triflate) ions act as spectators. The ytterbium cation is believed to be hydrated by the water present in the aqueous nitric acid solution. Nitric acid binds strongly to the hydrated ytterbium cation, as shown in Equation 1. A proton is generated, as shown in Equation 2, by the strong polarizing effect of the metal. Nitronium ion is then formed by the process shown in Equation 3. Although the nitronium ion may serve as the active electrophilic species, it is more likely that a nitronium carrier, such as the

¹ C. Braddock, Novel Recyclable Catalysts for Atom Economic Aromatic Nitration. *Green Chemistry*, 3 (2001): G26–G32.

intermediate formed in Equation 2, may serve as the electrophile. In any case, the reaction yields a nitrated aromatic compound.

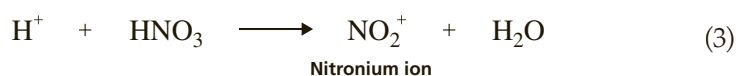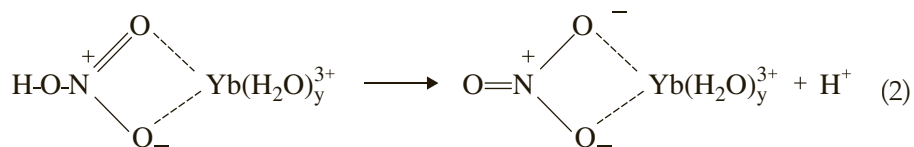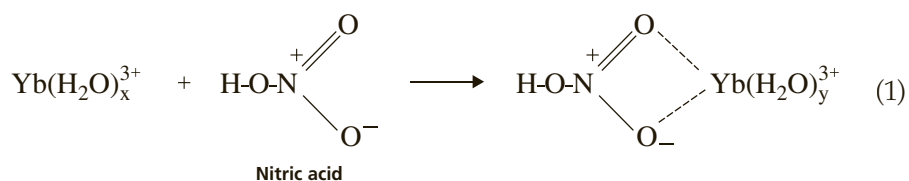

In this experiment, you will nitrate an aromatic substrate and analyze the composition of the mixture obtained by gas chromatography-mass spectrometry (GC-MS). In some cases, starting material will also be present in the mixture. You should be able to explain, mechanistically, why the observed products are obtained from the reaction.

REQUIRED READING

<i>Review:</i>	Technique 7	Reaction Methods, Sections 7.2 and 7.10
	Technique 7	Section 7.11 (optional)
	Technique 12	Extractions, Separations, and Drying Agents, Sections 12.4 and 12.9
	Technique 22	Gas Chromatography
<i>New:</i>	Essay	Green Chemistry
	Technique 28	Mass Spectrometry

SPECIAL INSTRUCTIONS

Some of the nitrated products may be toxic. All work should be conducted in a fume hood. Wear protective gloves to avoid skin contact with the nitrated products.

SUGGESTED WASTE DISPOSAL

The aqueous layer contains the catalyst, ytterbium triflate. Do not discard it. Instead, recycle the catalyst for future use by evaporating water on a hot plate. Transfer the colorless solid to a storage container or submit it to the instructor. If the material is highly colored, ask your instructor for advice. If the solvent, 1,2-dichloroethane, has been recovered using a rotary evaporator, pour it into a container so that it can be recycled.

NOTES TO THE INSTRUCTOR

It is suggested that each pair of students select a different substrate from the list provided. In most cases, the reaction will not go to completion, and as expected, will provide isomeric products. For example, toluene yields the expected *ortho* and *para* products, but a small amount of *meta* product is also formed. The products are analyzed by GC-MS. This experiment provides an excellent opportunity to discuss mass spectrometry because most of the compounds yield abundant molecular ions. The products are identified by searching the National Institute of Standards and Technology (NIST) database. Although it is best to search the database to identify the compounds, the experiment can also be conducted with gas chromatography. If this is done, one can usually assume that the nitro compounds will emerge in the following order: *ortho*, *meta*, and *para*. Adequate separations can be achieved on a GC-MS instrument using a J & W DB-5MS or Varian CP-Sil 5CB capillary column (30 m, 0.25-mm ID, 0.25 μ m). Set the injector temperature at 260°C. The column oven conditions are the following: start at 60°C (hold for 1 min), increase to 280°C at 20°C/min (12 min), and then hold at 280°C (4.5 min). Each run takes about 17 minutes. The helium flow rate is 1 mL/min. The mass range is set for 40 to 400 m/e.

PROCEDURE

Select one of the following aromatic substrates:

Toluene	Biphenyl
Butylbenzene	4-Methylbiphenyl
Isopropylbenzene	Diphenylmethane
<i>tert</i> -Butylbenzene	Phenylacetic acid
<i>ortho</i> -Xylene	Fluorobenzene
<i>meta</i> -Xylene	Iodobenzene
<i>para</i> -Xylene	Naphthalene
Anisole	Fluorene
1,2-Dimethoxybenzene (Veratrole)	Acetanilide
1,3-Dimethoxybenzene	Phenol
1,4-Dimethoxybenzene	α -Naphthol
4-Methoxytoluene	β -Naphthol

Place 0.375 g ytterbium (III) trifluoromethanesulfonate hydrate catalyst (ytterbium triflate) into a 25-mL round-bottom flask. Add 10 mL of 1,2-dichloroethane solvent followed by 0.400 mL of concentrated nitric acid (automatic pipette). Add two boiling stones to the flask. To this solution, weigh out and add approximately 6 millimoles of the aromatic substrate. Connect the round-bottom flask to a reflux condenser and clamp it into place on a ring stand. Use a very slow flow of water through the condenser. With a hot plate, heat the mixture to reflux for 1 hour.

After refluxing the mixture for 1 hour, allow the mixture to cool to room temperature and add 8 mL of water. Transfer the mixture into a separatory funnel. Shake the mixture gently and allow the two phases to separate. Drain the organic layer (bottom layer) into a 25-mL Erlenmeyer flask. Dry the organic layer with a small scoop of anhydrous magnesium sulfate (about 0.5 g). If a rotary evaporator is available, transfer the organic layer to a preweighed 50-mL round-bottom flask for removal of solvent. The apparatus allows the possibility of recovering most of the 1,2-dichloroethane. When the solvent has been removed, remove the flask and weigh it.

Alternatively, the solvent may be removed by using the apparatus shown in Technique 7, Figure 7.17C. Transfer the dried organic layer to a preweighed 125-mL filter flask. Add a melting-point capillary tube to the flask (open end down) and then cork the top. The melting-point capillary tube helps speed the evaporation process. Connect the sidearm of the filter flask to the house vacuum system or aspirator, using a trap cooled in ice. There will be a cooling effect while the evaporation takes place, so you will need to heat the flask gently (lowest setting on a hot plate). Most of your solvent should be evaporated in less than 1 hour, under vacuum and with gentle heating. Weigh the filter flask.

The aqueous layer remaining in the separatory funnel contains the ytterbium catalyst. Pour the aqueous layer from the top of the separatory funnel into a preweighed 50-mL Erlenmeyer flask. Completely evaporate the water on a hot plate. Weigh the flask to determine how much catalyst you were able to recover. Place the catalyst in a container that holds the recycled catalyst that will be reused in other classes.

Unless instructed otherwise, prepare a sample for analysis by GC-MS by dissolving two drops of the mixture of nitrated aromatic compounds in about 1 mL of methylene chloride. These samples will be run using automation software on the GC-MS system.

When your sample has been run, you will have an opportunity to search the NIST mass spectral library to determine the structures of the product(s) of the nitration. Determine the structures of the product(s) and the percentages of each component. There will likely be starting material left in the reaction mixture. It would be of interest to see how your product ratios compare to the values obtained from the literature (see References).

REFERENCES

- Braddock, C. Novel Recyclable Catalysts for Atom Economic Aromatic Nitration. *Green Chem.* **2001**, 3, G26–G32.
- Schofield, K. *Aromatic Nitration*; Cambridge University Press: London, 1980.
- Waller, F. J.; Barrett, G. M.; Braddock, D. C.; Ramprasad, D. Lanthanide (III) Triflates as Recyclable Catalysts for Atom Economic Aromatic Nitration. *Chem. Comm.* **1997**, 613–614.

QUESTIONS

1. Interpret the mass spectrum of the compounds formed in the nitration of your aromatic substrate.
2. Draw a mechanism that explains how the nitro-substituted aromatic products observed in your reaction were formed.


Reduction of Ketones Using Carrots as Biological Reducing Agents

Green chemistry

Use of a biological reducing agent

Reduction of a ketone to an alcohol


A very common reaction in organic chemistry is the reduction of a ketone to a secondary alcohol.

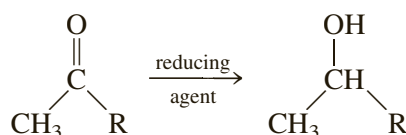

The most widely used reducing agents include lithium aluminum hydride, sodium borohydride (see Experiment 33), and catalytic hydrogenation. The reaction takes place in an organic solvent, such as diethyl ether or methanol.


Biological reducing agents can also be used to bring about the reduction of a ketone to a secondary alcohol. The reduction of the carbonyl group of ethyl acetoacetate (Experiment 29) is carried out using baker's yeast as a reducing medium. In this experiment, grated carrot is used to bring about a similar reaction. This type of reaction is an example of a green chemistry application, because water is the only solvent and the principal reagent is a common garden vegetable.

In each of these biological reduction experiments, an organic molecule is used by the biological system as the actual reducing agent. This reducing agent is **nicotinamide adenine dinucleotide (NADH)**. NADH acts as a *cofactor*; its chemical properties are expressed in coordination with an enzyme, which regulates the process.

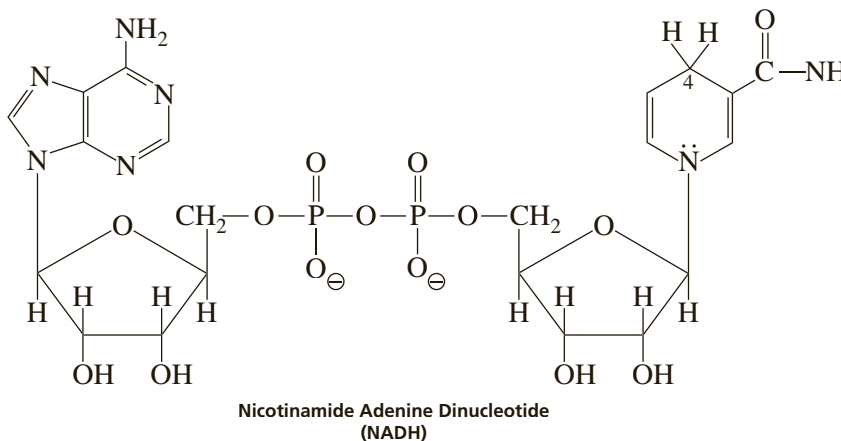


While the structure of NADH may seem overly complex, it is only necessary to focus on the nicotinamide ring—specifically on the hydrogen atoms attached to C4. This is the actual reactive site of the NADH molecule; the rest of the structure is important for enzyme–substrate binding, water solubility, ease of transport through cell walls, and so on.

In the biological reduction of a ketone, one of the hydrogens at C4 of the nicotinamide ring is transferred *with its pair of electrons*, in the form of a hydride, to the carbonyl carbon of the ketone. Note that the hydride is acting as a *nucleophile* as it attacks the carbonyl carbon.

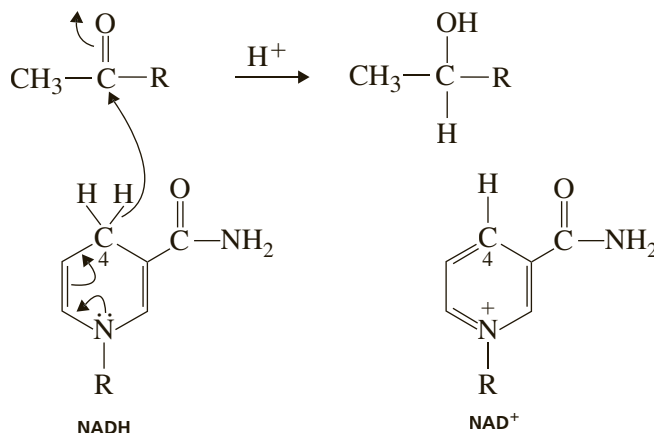

In the process of reducing the ketone, NADH is oxidized to NAD⁺. This reaction is energetically favorable, because the aromatic property of the pyridine ring is restored—a gain in *resonance energy*.

In this experiment, the biological source of NADH will be a common garden carrot. The reaction is:

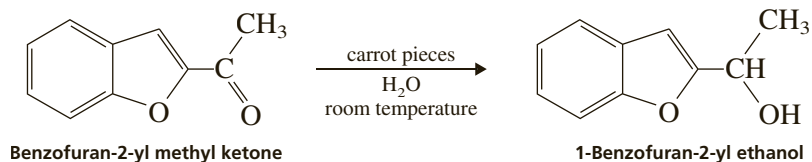

The results of the reduction will be analyzed by infrared spectroscopy. While we might expect this reduction to be stereoselective, the scale of the reaction used here does not permit an optical purity analysis by polarimetry.

REQUIRED READING

<i>Review:</i> Essay	Green Chemistry
Technique 8	Filtration
Technique 12	Extractions, Separations, and Drying Agents
Technique 25	Infrared Spectroscopy

SPECIAL INSTRUCTIONS

This experiment must be allowed to stand aside for a period of at least 24 hours. Another experiment can be conveniently co-scheduled with this one.

SUGGESTED WASTE DISPOSAL

The carrot residue may be safely disposed of in the trash. Diethyl ether solvent can be recovered after the evaporation step if a rotary evaporator is available.

PROCEDURE

Grate one fresh carrot to obtain approximately 25 g of grated carrot. Wash this carrot material with distilled water. Weigh the grated carrot in a 150-mL Erlenmeyer flask, and add 75 mL of distilled water and a magnetic stirring bar. Add 50 mg of benzofuran-2-yl methyl ketone to the flask, stopper it with a cork, and clamp it in position above a magnetic stirrer. Allow the mixture to stir for at least 24 hours. Be sure to clamp the flask so that there is some space between the bottom of the flask and top of the magnetic stirrer. This is to avoid any heating from the stirrer motor, which may stop the reaction.

After the stirring has been stopped, filter the reaction mixture through a layer of cheesecloth to remove the larger chunks of carrot. Remove the remaining carrot residue by vacuum filtration using a Hirsch funnel (see Technique 8, Section 8.3).

Extract the filtrate three times with 10-mL portions of diethyl ether. Dry the ether extract over anhydrous magnesium sulfate. Transfer the dried solution into a clean flask, and remove the ether solvent by evaporation (if a rotary evaporator is available, it is best to use it).

Determine the infrared spectrum of the product as a neat liquid (see Technique 25, Section 25.2). You should be able to observe the extent of the reduction by noting the disappearance of the carbonyl stretching peak at about 1700 cm^{-1} and the appearance of a strong O—H stretching peak at about 3450 cm^{-1} . Be sure to submit your spectra with your laboratory report.

REFERENCE

Ravia, S.; Gamemara, D.; Schapiro, V.; Bellomo, A.; Adum, J.; Seoane, G.; Gonzalez, D. Enantioselective Reduction by Crude Plant Parts: Reduction of Benzofuran-2-yl Methyl Ketone with Carrot (*Daucus carota*) Bits. *J. Chem. Educ.* **2006**, *83*, 1049–1051.

Resolution of (\pm)- α -Phenylethylamine and Determination of Optical Purity

Resolution of enantiomers

Use of a separatory funnel

Polarimetry

Chiral gas chromatography

NMR spectroscopy


Chiral resolving agent

Diastereomeric methyl groups

Although racemic (\pm)- α -phenylethylamine is readily available from commercial sources, the pure enantiomers are more difficult to obtain. In this experiment, you will isolate one of the enantiomers, the levorotatory one, in a high state of optical purity (large enantiomeric excess). A **resolution**, or separation, of enantiomers will be performed, using (+)-tartaric acid as the resolving agent.

Resolution of Enantiomers

The resolving agent to be used is (+)-tartaric acid, which forms diastereomeric salts with racemic α -phenylethylamine. The important reactions for this experiment follow.

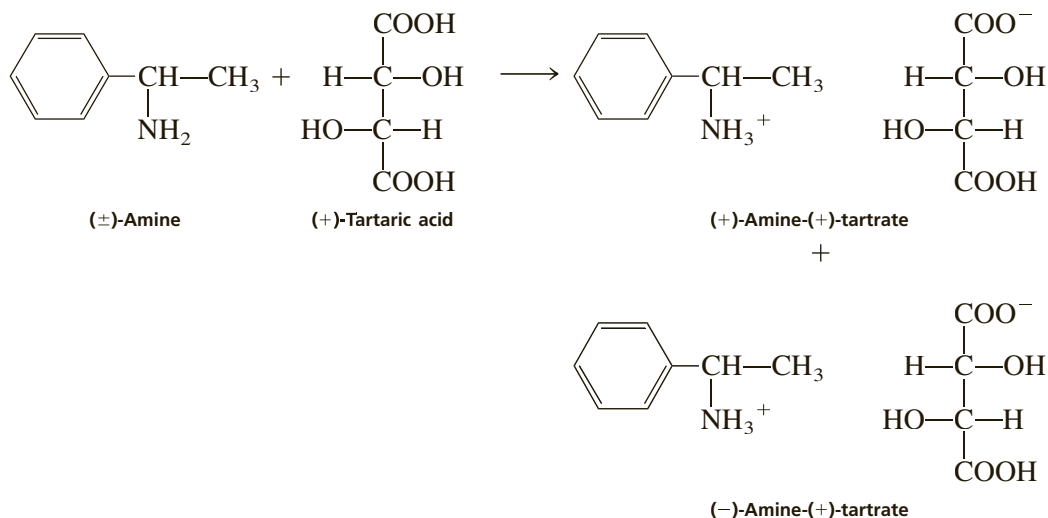

Optically pure (+)-tartaric acid is abundant in nature. It is frequently obtained as a by-product of winemaking. The separation depends on the fact that diastereomers usually have different physical and chemical properties. The (–)-amine-(+)–tartrate salt has a lower solubility than its diastereomeric counterpart, the (+)-amine-(+)–tartrate salt. With some care, the (–)-amine-(+)–tartrate


salt can be induced to crystallize, leaving (+)-amine-(+)-tartrate in solution. The crystals are removed by filtration and purified. The (–)-amine can be obtained from the crystals by treating them with base. This breaks apart the salt by removing the proton, and it regenerates the free, unprotonated (–)-amine.

A polarimeter will be used to measure the observed rotation, α , of the resolved amine sample. From this value, you will calculate the specific rotation $[\alpha]_D$ and the optical purity (enantiomeric excess) of the amine. You will then calculate the percentages of each of the enantiomers present in the resolved sample. The (*S*)- α -phenylethylamine predominates in the sample. An optional chiral gas chromatographic method may be used to directly determine the percentages of each of the enantiomers in the sample.


NMR Determination of Optical Purity

An alternate means of determining the optical purity of the sample makes use of NMR spectroscopy (see Experiment 32B). A group attached to a stereogenic (chiral) carbon normally has the same chemical shift whether that carbon has either an *R* or *S* configuration. However, that group can be made diastereomeric in the NMR spectrum (have different chemical shifts) when the racemic parent compound is treated with an optically pure chiral resolving agent to produce diastereomers. In this case, the group is no longer found in two enantiomers but, rather, in two different diastereomers, and its chemical shift will be different in each environment.

In this experiment, the partly resolved amine (containing both *R* and *S* enantiomers) is mixed with optically pure (*S*)-(+)-*O*-acetylmandelic acid in an NMR tube containing CDCl_3 . Two diastereomers are formed:

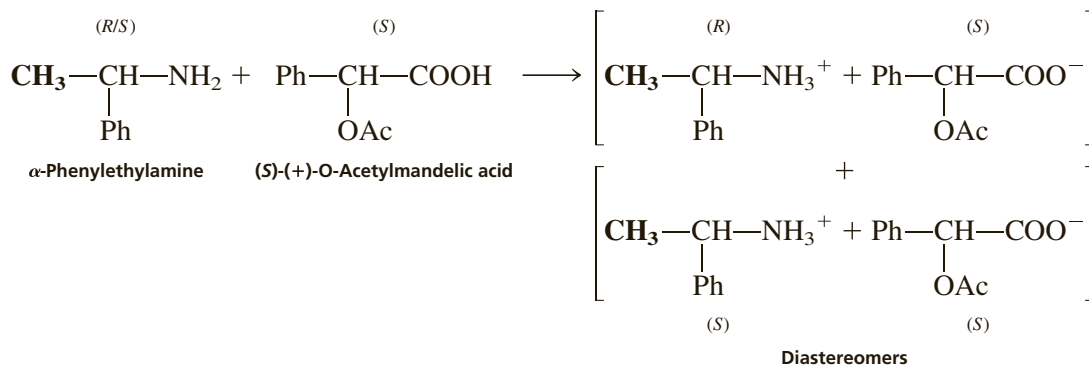

The methyl groups in the amine portions of the two diastereomeric salts are attached to a stereocenter, (*S*) in one case and (*R*) in the other. As a result, the methyl groups themselves become diastereomeric, and they have different chemical shifts. In this case, the (*R*) isomer is downfield, and the (*S*) isomer is upfield. These methyl groups appear at approximately (varies) 1.1 and 1.2 ppm, respectively, in the proton NMR spectrum of the mixture. Because the methyl groups are adjacent to a methine (CH) group, they appear as doublets. These doublets may be integrated in order to determine the percentage of the (*R*) and (*S*) amines in the resolved α -phenylethylamine. In the example, the NMR spectrum was determined with a mixture made by dissolving equal quantities (50:50 mixture) of the original unresolved (\pm)- α -phenylethylamine and a student's resolved product, which contained predominantly (*S*)-(+)- α -phenylethylamine.

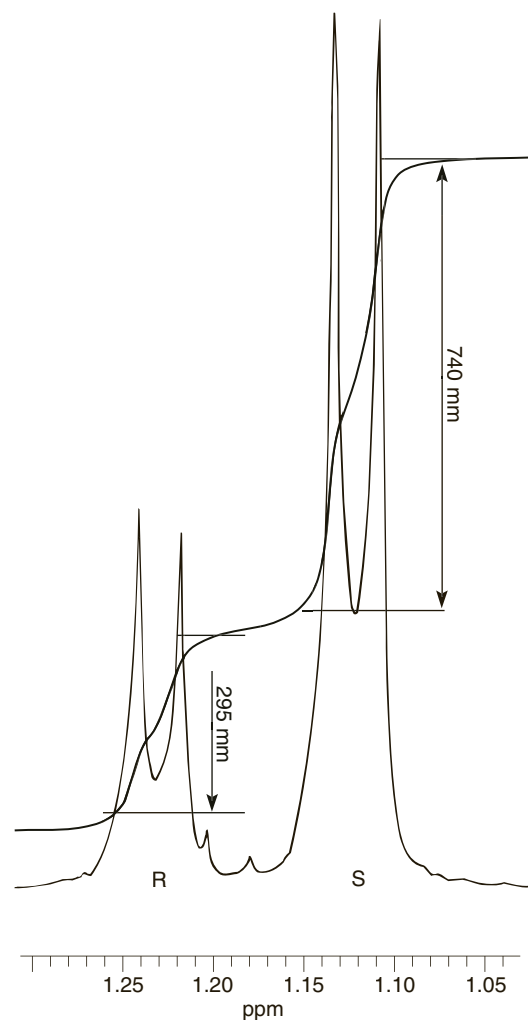

300-MHz Spectrum of a 50:50 mixture of resolved and unresolved α -phenylethylamine in CDCl_3 . The chiral resolving agent (S)-(+)-O-acetylmandelic acid was added.

32A EXPERIMENT 32A

Resolution of (\pm)- α -Phenylethylamine

In this procedure, you will resolve racemic (\pm)- α -phenylethylamine, using (+)-tartaric acid as the resolving agent.

REQUIRED READING

Review: Technique 8 Section 8.3
Technique 12 Sections 12.4, 12.8, 12.9
Technique 23
Technique 22 (optional)

SPECIAL INSTRUCTIONS

α -Phenylethylamine readily reacts with carbon dioxide in the air to form a white solid, the *N*-carboxyl amine derivative. Every effort should be taken to avoid prolonged exposure of the amine to air. Be sure to close the bottle tightly after you have measured the rotation of your amine and be sure to place your sample quickly into the flask where you will perform the resolution. This flask should also be stoppered. Use a cork stopper because a rubber stopper will dissolve somewhat and discolor your solution. The crystalline salt will not react with carbon dioxide until you decompose it to recover the resolved amine. Then, you must be careful once again.

The observed rotation for a sample isolated by a single student may be only a few degrees, which limits the precision of the optical purity determination. Better results can be obtained if four students combine their resolved amine products for the polarimetric analysis. If you have allowed your amine to have excessive exposure to air, the polarimetry solution may be cloudy. This will make it difficult to obtain an accurate determination of the optical rotation.

SUGGESTED WASTE DISPOSAL

Place the mother liquor solution from the crystallization, which contains (+)- α -phenylethylamine, (+)-tartaric acid, and methanol, in the special container provided for this purpose. Aqueous extracts will contain tartaric acid, dilute base, and water; they should be placed in the container designated for aqueous wastes. When you are finished with polarimetry, depending on the wishes of your instructor, you should either place your resolved (S)-(-)- α -phenylethylamine in a special container marked for this purpose or you should submit it to your instructor in a suitably labeled container that includes the names of those people who have combined their samples.

PROCEDURE

NOTE TO THE INSTRUCTOR

This experiment is designed for students to work individually, but to combine their products with three other students for polarimetry.

Preparations Place 7.8 g of L-(+)-tartaric acid and 125 mL of methanol in a 250-mL Erlenmeyer flask. Heat this mixture on a hot plate until the solution is nearly boiling. Slowly add 6.25 g of racemic α -phenylethylamine (α -methylbenzylamine) to this hot solution.

CAUTION



At this step, the mixture is likely to froth and boil over.

Crystallization

Stopper the flask and let it stand overnight. The crystals that form should be prismatic. If needles form, they are not optically pure enough to give a complete resolution of the enantiomers; *prisms must form*. Needles should be dissolved (by

careful heating) and cooled slowly to crystallize once again. When you recrystallize, you can “seed” the mixture with a prismatic crystal, if one is available. If it appears that you have prisms but that they are overgrown (covered) with needles, the mixture may be heated. The mixture may be heated until *most* of the solid has dissolved. The needle crystals dissolve easily, and usually a small amount of the prismatic crystals remains to seed the solution. After dissolving the needles, allow the solution to cool slowly and form prismatic crystals from the seeds.

Workup

Filter the crystals, using a Büchner funnel (see Technique 8, Section 8.3, and Figure 8.5), and rinse them with a few portions of cold methanol. Partially dissolve the crystalline amine-tartrate salt in 25 mL of water, add 4 mL of 50% sodium hydroxide, and extract this mixture with three 10-mL portions of methylene chloride using a separatory funnel (see Technique 12, Section 12.4). Combine the organic layers from each extraction in a stoppered flask and dry them over about 1 g of anhydrous sodium sulfate for about 10 minutes.

Two different methods should be considered for removing the solvent. Ask your instructor which method you should use. **Method 1** involves using a rotary evaporator to remove the solvent. If you are employing this method, preweigh a 100-mL *round*-bottom flask, and decant the methylene chloride solution containing the amine into the flask. Ask your instructor to demonstrate the use of the rotary evaporator. A liquid remains after the solvent has been removed. You may need to increase the temperature of the water bath to ensure that all of the solvent has been removed. About 2 or 3 mL of the liquid amine should remain. Proceed to the Yield Calculation and Storage section below.

If your instructor asks you to use **Method 2**, proceed as follows. While the solution is drying over anhydrous sodium sulfate, preweigh a clean, dry 50-mL Erlenmeyer flask. Decant the dried solution into the flask and evaporate the methylene chloride on a hot plate (about 60°C) in a hood. A stream of nitrogen or air should be directed into the flask to increase the rate of evaporation. When the volume of liquid reaches about 2 or 3 mL total, you should carefully insert a hose attached to the house vacuum or aspirator system to remove any remaining methylene chloride. The hose should be inserted into the neck of the flask. Note that the desired product is a **liquid**. Some solid amine carbonate may start to form on the sides of the flask during the course of the evaporation. This undesired solid is more likely to form if you prolong the heating operation. You will want to take care to avoid the formation of this white solid if at all possible. If you do obtain a cloudy solution or solids are present, transfer the material to a centrifuge tube and centrifuge the sample. Then remove the clear liquid for the polarimetry part of this experiment.

Yield Calculation and Storage

Stopper the flask and weigh it to determine the yield. Also calculate the percentage yield of the (S)-(-)-amine based on the amount of the racemic amine you started with.

Polarimetry

Combine your product with the products obtained by three other students. If anyone's product is highly colored or if a large amount of solid is present, do not use it. If the amine is a little cloudy or if there is just a small amount of solid present, transfer the sample to a small centrifuge tube (microcentrifuge tubes work well here) and centrifuge the sample for about 5 minutes. Remove the *clear* liquid with

a Pasteur pipette to avoid drawing up any solid into the pipette and fill a pre-weighed 10-mL volumetric flask. You will not get good results with the polarimeter if the amine is cloudy or if there are suspended solids present in your amine, so be careful to avoid transferring any solid.

Weigh the flask to determine the weight of amine and calculate the density (concentration) in grams per milliliter. You should obtain a value of about 0.94 g/mL. This should give you a sufficient amount of material to proceed with the polarimetry measurements that follow without diluting your sample. If, however, your combined products do not amount to more than 10 mL of the amine, you may have to dilute your sample with methanol (check with your instructor).

If you have less than 10 mL of product, weigh the flask to determine the amount of the amine present. Then fill the volumetric flask to the mark with absolute methanol and mix the solution thoroughly by inverting 10 times. The concentration of your solution in grams per milliliter is easily calculated.

Transfer the solution to a 0.5-dm polarimeter tube and determine its observed rotation. Your instructor will show you how to use the polarimeter. Report the values of the observed rotation, specific rotation, and optical purity (enantiomeric excess) to the instructor. The published value for the specific rotation is $[\alpha_D^{22}] = -40.3^\circ$. Calculate the percentage of *each* of the enantiomers in the sample (see Technique 23, Section 23.5), and include the figures in your report.

Due to the presence of some methylene chloride in the sample of the chiral amine, you may obtain low rotation values from polarimetry. Because of this, your calculated value of the optical purity (enantiomeric excess) and percentages of the enantiomers will be in error. The percentages of the enantiomers obtained from the optional chiral gas chromatography experiment below should provide more accurate percentages of each of the stereoisomers.

Chiral Gas Chromatography (Optional)

Chiral gas chromatography will provide a direct measure of the amounts of each stereoisomer present in your resolved α -phenylethylamine sample. A Varian CP-3800 equipped with J & P (Agilent) Cyclosil B capillary column (30 m, 0.25-mm ID, 0.25 μ m) provides an excellent separation of (*R*) and (*S*)-enantiomers.

Set the FID detector at 270°C and the injector temperature at 250°C. The initial split ratio should be set at 150:1 and then changed to 10:1 after 1.5 minutes. Set the oven temperature at 100°C and hold at that temperature for 25 minutes. The helium flow rate is 1 mL/min. The compounds elute in the following order: (*R*)- α -phenylethylamine (17.5 min) and (*S*)-enantiomer (18.1 min). Your observed retention times may vary from those given here, but the order of elution will be the same. Because the peaks overlap slightly, you may not observe a distinct peak for the (*R*)-enantiomer. Instead, you may observe a shoulder for the (*R*)-enantiomer peak on the side of the large peak for the (*S*)-enantiomer. If you are able to see the (*R*)-enantiomer, integrate the area under the peaks to obtain the percentages of each of the enantiomers in your sample and compare your results to those obtained with the polarimeter. It should be noted that the resolution process used in this experiment is highly selective for the (*S*)-enantiomer. That is the good news; the bad news is that you may have such a pure (*S*)- α -phenylethylamine sample that you will not be able to obtain percentages from the analysis on the chiral column.

32B EXPERIMENT 32B

Determination of Optical Purity Using NMR and a Chiral Resolving Agent

In this procedure, you will use NMR spectroscopy with the chiral resolving agent (S)-(+)-O-acetylmandelic acid to determine the optical purity of the (S)-(-)- α -phenylethylamine you isolated in Experiment 32A.

REQUIRED READING

New: Technique 26 Nuclear Magnetic Resonance Spectroscopy

SPECIAL INSTRUCTIONS

Be sure to use a clean Pasteur pipette whenever you remove CDCl_3 from its supply bottle. Avoid contaminating the stock of NMR solvent. Also be sure to fill and empty the pipette several times before attempting to remove the solvent from the bottle. If you bypass this equilibration technique, the volatile solvent may squirt out of the pipette before you can transfer it successfully to another container.

SUGGESTED WASTE DISPOSAL

When you dispose of your NMR sample, which contains CDCl_3 , place it in the container designated for halogenated wastes.

PROCEDURE

Using a small test tube, weigh approximately 0.05 mmole (0.006 g, $MW = 121$) of your resolved amine by adding it from a Pasteur pipette. Cork the test tube to protect it from atmospheric carbon dioxide. Carbon dioxide reacts with the amine to form an amine carbonate (white solid). Using a weighing paper, weigh approximately 0.06 mmole (0.012 g, $MW = 194$) of (S)-(+)-O-acetylmandelic acid and add it to the amine in the test tube. Using a clean Pasteur pipette, add about 0.25 mL of CDCl_3 to dissolve everything. If the solid does not completely dissolve, you can mix the solution by drawing it several times into your Pasteur pipette and redelivering it back into the test tube. When everything is dissolved, transfer the mixture to an NMR tube using a Pasteur pipette. Using a clean Pasteur pipette, add enough CDCl_3 to bring the total height of the solution in the NMR tube to 50 mm.

Determine the proton NMR spectrum, preferably at 300 MHz, using a method that expands and integrates the peaks of interest. Using the integrals, calculate the



percentages of the *R* and *S* isomers in the sample and its optical purity.¹ Compare your results from this NMR determination to those you obtained by polarimetry (Experiment 32A).

REFERENCES

- Ault, A. Resolution of D, L- α -Phenylethylamine. *J. Chem. Educ.* **1965**, *42*, 269.
 Jacobus, J.; Raban, M. An NMR Determination of Optical Purity. *J. Chem. Educ.* **1969**, *46*, 351.
 Parker, D.; Taylor, R. J. Direct ¹H NMR Assay of the Enantiomeric Composition of Amines and β -Amino Alcohols Using O-Acetyl mandelic Acid as a Chiral Solvating Agent. *Tetrahedron* **1987**, *43* (22), 5451.

QUESTIONS

- Using a reference textbook, find examples of reagents used in performing chemical resolutions of acidic, basic, and neutral racemic compounds.
- Propose methods of resolving each of the following racemic compounds

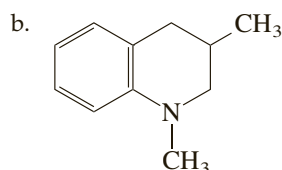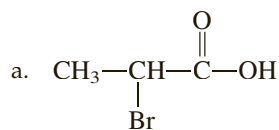

- Explain how you would proceed to isolate (*R*)-(+)- α -phenylethylamine from the *mother liquor* that remained after you crystallized (*S*)-(–)- α -phenylethylamine.
- What is the white solid that forms when α -phenylethylamine comes in contact with carbon dioxide? Write an equation for its formulation.
- Which method, polarimetry or NMR spectroscopy, gives the more accurate results in this experiment? Explain.
- Draw the three-dimensional structure of (*S*)-(–)- α -phenylethylamine.
- Draw the three-dimensional structure of the diastereomer formed when (*S*)-(–)- α -phenylethylamine is reacted with (*S*)-(+)-O-acetylmandelic acid.

¹ *Note to the Instructor:* In some cases, the resolution is so successful that it is very difficult to detect the doublet arising from the (*R*)-(+)- α -phenylethylamine + (*S*)-(+)-O-acetylmandelic acid diastereomer. If this occurs, it is useful to have the students add a single drop of racemic α -phenylethylamine to the NMR tube and re-determine the spectrum. In this way, both diastereomers can be clearly seen.

An Oxidation–Reduction Scheme: Borneol, Camphor, Isoborneol¹

Green chemistry

Oxone[®] oxidation

Sodium borohydride reduction

Sodium borohydride reduction


Sublimation (optional)

Stereochemistry

Gas chromatography

Spectroscopy (infrared, proton NMR, carbon-13 NMR)

Computational chemistry (optional)

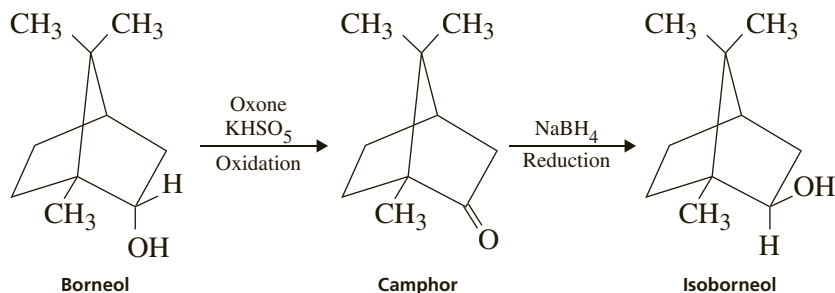

This experiment will illustrate the use of a “green” oxidizing agent, Oxone[®], for converting a secondary alcohol (borneol) to a ketone (camphor). The camphor is then reduced by sodium borohydride to give the *isomeric* alcohol, isoborneol. The spectra of borneol, camphor, and isoborneol will be compared to detect structural differences and to determine the extent to which the final step produces an alcohol isomeric with the starting material, borneol.


OXIDATION OF BORNEOL WITH OXONE[®]

The oxidizing agent in this reaction is Oxone[®], which is a “triple salt” (2KHSO₅·KHSO₄·K₂SO₄). Oxone[®] is manufactured by DuPont and is sometimes used in swimming pools in conjunction with chlorine. Oxone[®] reacts with organic contaminants that are introduced by swimmers and the environment. It can also be used in the laboratory as the oxidizing agent in oxidation reactions of organic compounds.

¹This experiment was adapted with permission from “Oxidation of Borneol to Camphor Using Oxone[®] and Catalytic Sodium Chloride: A Green Experiment for the Undergraduate of Organic Chemistry Laboratory” by Patrick T. Lang; Andrew M. Harned; Jane E. Wissinger. *J. Chem. Educ.* **2011**, *88*, 652–656.

It is a much less toxic alternative to the use of chromium reagents, which have traditionally been used in these oxidation reactions. The active oxidizing agent in Oxone® is potassium peroxydisulfate, $\text{K}_2\text{S}_2\text{O}_8$. Chloride ions, which are added to the reaction mixture as sodium chloride, are oxidized by potassium peroxydisulfate to molecular chlorine. In aqueous solution, some chlorine is converted to hypochlorous acid (HOCl). The actual oxidizing agent in the conversion of borneol to camphor is most likely hypochlorous acid.


Although the mechanism is not fully understood, there is evidence that an alkyl hypochlorite intermediate is produced, which then gives the product via an E_2 elimination:

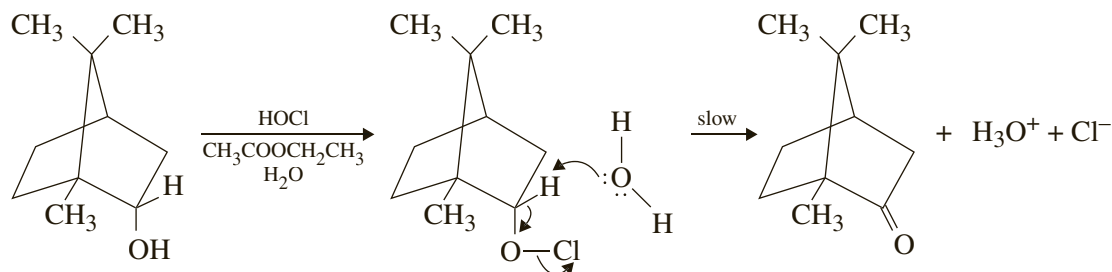

REDUCTION OF CAMPHOR WITH SODIUM BOROHYDRIDE

Metal hydrides (sources of H^-) of the Group III elements, such as lithium aluminum hydride, LiAlH_4 , and sodium borohydride, NaBH_4 , are widely used in reducing carbonyl groups. Lithium aluminum hydride, for example, reduces many compounds containing carbonyl groups such as aldehydes, ketones, carboxylic acids, esters, or amides, whereas sodium borohydride reduces only aldehydes and ketones. The reduced reactivity of borohydride allows it to be used even in alcohol and water solvents, whereas lithium aluminum hydride reacts violently with these solvents to produce hydrogen gas and thus must be used in nonhydroxylic solvents. In the present experiment, sodium borohydride is used because it is easily handled, and the results of reductions using either of the two reagents are essentially the same. The same care required for lithium aluminum hydride in keeping it away from water need not be taken for sodium borohydride.

The mechanism of action of sodium borohydride in reducing a ketone is as follows:

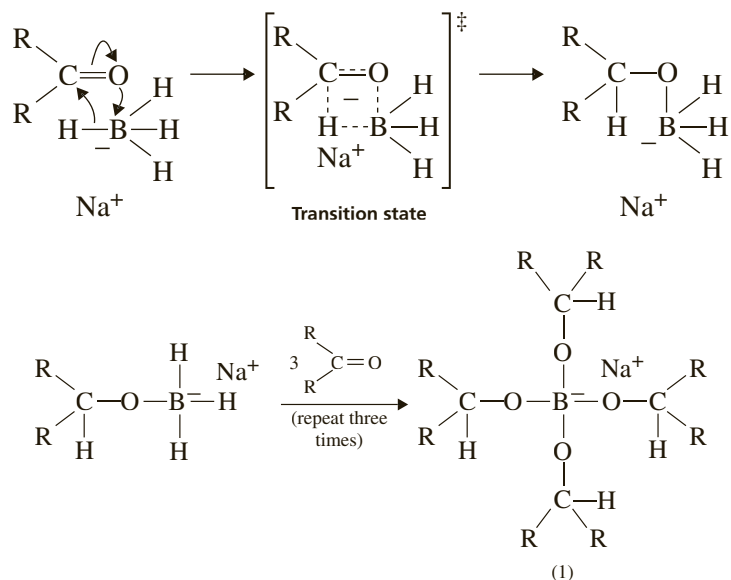


Note in this mechanism that all four hydrogen atoms are available as hydrides (H^-), and thus one mole of borohydride can reduce four moles of ketones. All of the steps are irreversible. Usually, excess borohydride is used because there is uncertainty regarding its purity and because some of it reacts with the solvent.

Once the final tetraalkoxyboron compound (1) is produced, it can be decomposed (along with excess borohydride) at elevated temperatures as shown:




The stereochemistry of the reduction is very interesting. The hydride can approach the camphor molecule more easily from the bottom side (**endo** approach) than from the top side (**exo** approach). If attack occurs at the top, a large steric repulsion is created by one of the two **geminal** methyl groups. Geminal methyl groups are groups that are attached to the same carbon. Attack at the bottom avoids this steric interaction.

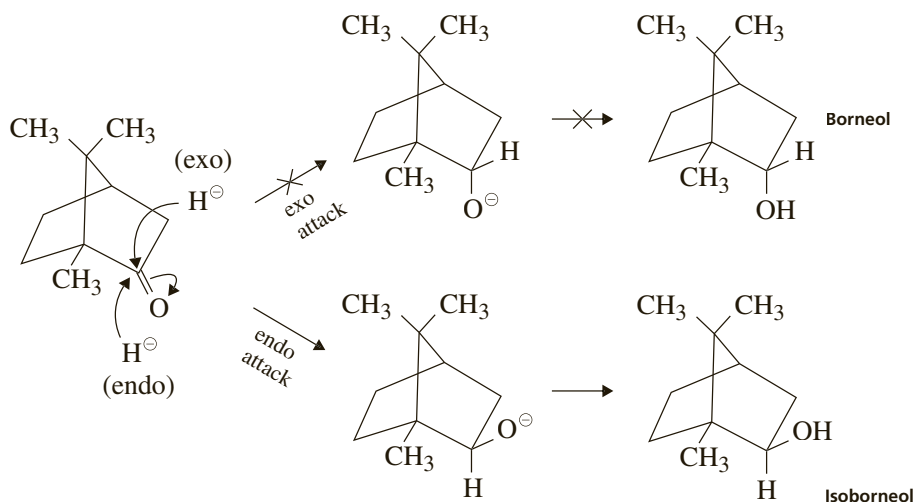


It is expected, therefore, that **isoborneol**, alcohol produced from the attack at the *least*-hindered position, will *predominate but will not be the exclusive product* in the final reaction mixture. The percentage composition of the mixture can be determined by spectroscopy.

It is interesting to note that when the methyl groups are removed (as in 2-norbornanone), the top side (**exo** approach) is favored, and the opposite stereochemical result is obtained. Again, the reaction does not give exclusively one product.

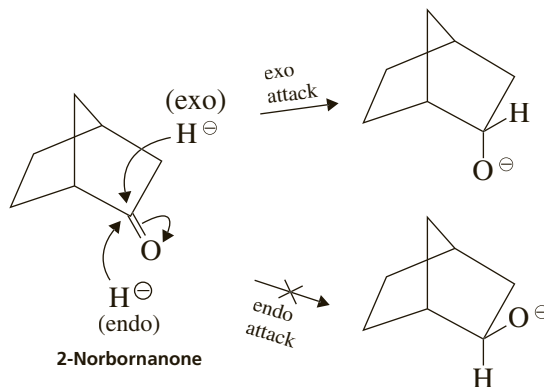

	NaBH_4	LiAlH_4
<i>exo</i>	86%	89%
<i>endo</i>	14%	11%

Bicyclic systems such as camphor and 2-norbornanone react predictably according to steric influences. This effect has been termed **steric approach control**. In the reduction of simple acyclic and monocyclic ketones, however, the reaction seems to be influenced primarily by thermodynamic factors. This effect has been termed **product development control**. In the reduction of 4-*t*-butylcyclohexanone, the thermodynamically more stable product, is produced by product development control.

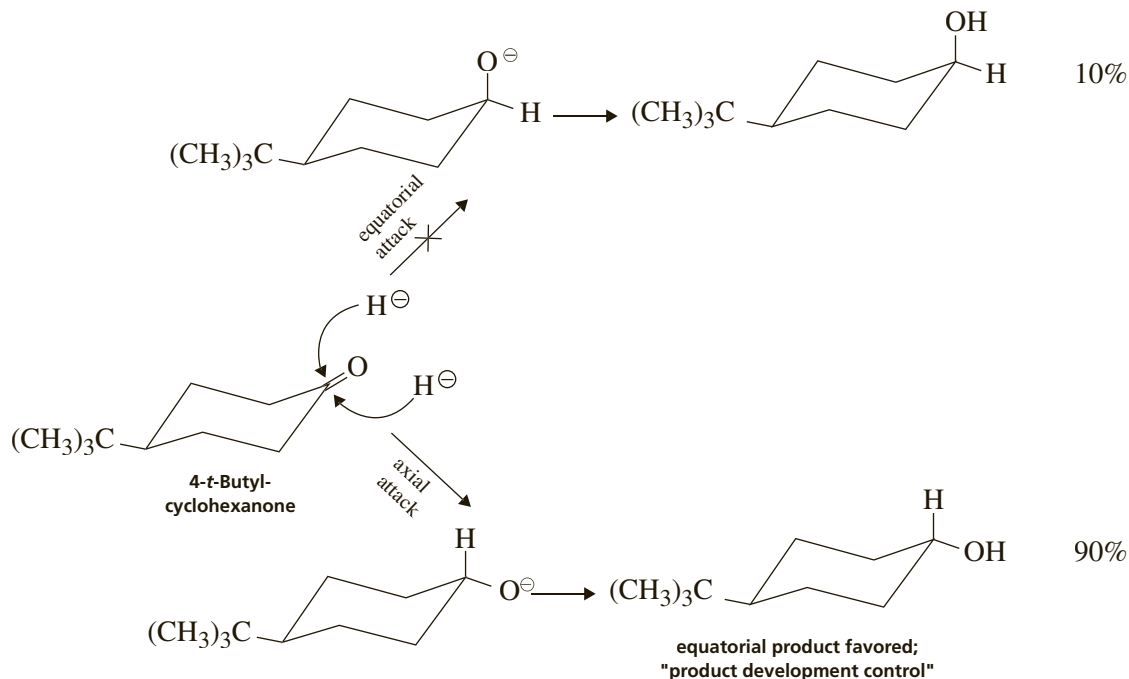

REQUIRED READING

<i>Review:</i>	Technique 6	Heating and Cooling Methods, Sections 6.1–6.3
	Technique 7	Reaction Methods, Sections 7.1–7.4 and 7.10
	Technique 8	Filtration, Section 8.3
	Technique 9	Physical Constants of Solids: The Melting Point, Sections 9.7 and 9.8
	Technique 12	Extractions, Separations, and Drying Agents Section 12.4
	Techniques 20, 22, 25, 26, and 27	
<i>New:</i>	Technique 17	Sublimation
	Essay	Green Chemistry
	Essay and Experiment 20	Computational Chemistry (optional)

SPECIAL INSTRUCTIONS

Oxone[®] is a strong oxidant. You should not inhale the dust. The reactants and products are all highly volatile and must be stored in tightly closed containers.

SUGGESTED WASTE DISPOSAL

The aqueous solutions obtained from the extraction steps should be placed in the aqueous waste container. Any leftover methanol may be placed in the nonhalogenated organic waste container.

NOTES TO THE INSTRUCTOR

We use Oxone[®] from Alfa Aesar, catalog # 89892.

The traditional test for the activity of sodium borohydride is to place a small amount of powdered material in some methanol, and heat it gently. The solution should bubble vigorously if the hydride is active. However, we have found this test to give unreliable results. The best way to test the activity of sodium borohydride is to run the reduction reaction with camphor and analyze the product with infrared spectroscopy.

Percentages of borneol and isoborneol can be determined by gas chromatography. Any gas chromatograph should be suitable for this determination. For example, a Gow-Mac 69-930 instrument with an 8-ft column of 10% Carbowax 20M, at 180°C and with a 40 mL/min helium flow rate, will give a suitable separation. The compounds elute in the following order: camphor (8 min), isoborneol (10 min), and borneol (11 min). A Varian CP-3800 with autosampler equipped with a J & W DB-5 or Varian CP-Sil 5CB capillary column (30 m, 0.25-mm ID, 0.25 μm) also provides a good separation. Set the injector temperature at 250°C. The column oven conditions are the following: start at 75°C (hold for 10 min), increase to 200°C at 35°C min, and then hold at 200°C (1 min). Each run takes about 15 minutes. The helium flow rate is 1 mL/min. The compounds elute in the following order: camphor (12.9 min), isoborneol (13.1 min), and borneol (13.2 min). An optional procedure is provided that involves computational chemistry.

PROCEDURE

Part A. Oxidation of Borneol to Camphor².

Setup of the Reaction. Clamp a 10-mL round-bottom flask to a ring stand with the flask sitting on a magnetic stirrer. Add a magnetic stir bar to the flask along with 0.5 g of borneol. After adding 2.0 mL of ethyl acetate to the flask, stir to dissolve the borneol. With continued stirring, add 1.2 g of Oxone[®], 0.04 g of sodium chloride, and 1.0 mL of deionized water. Stir the reaction mixture at room temperature for 50 minutes. Then add an additional 0.015 g of NaCl and stir for 10 more minutes. Note any changes in color during the entire procedure.

At this point the oxidation should be complete and the excess oxidant should be destroyed. Add 7.5 mL of deionized water to the flask and stir to dissolve most of the salts. Add about 0.1 g of solid sodium bisulfite to reduce the oxidants that remain. Test the aqueous layer for excess oxidant by removing a couple drops of the **bottom** layer using a Pasteur pipette. Put a drop of this liquid on a piece of starch-iodide paper. A blue-black color (positive test) indicates the presence of excess oxidant; add small amounts of sodium bisulfite if the aqueous layer tests positive, until a negative test is achieved (no color change).

Extraction of Camphor. Using a Pasteur pipette, transfer the contents to a centrifuge tube. Rinse the round-bottom flask with about 1 mL of ethyl acetate and transfer the remaining liquid to the centrifuge tube. Cap the tube and shake gently for 1 minute. Transfer the aqueous layer to a second centrifuge tube. (Note: It may be difficult to see the interface between the two layers, but it will be there!) Extract the aqueous layer twice more, with 2 mL of ethyl acetate each time. Combine the three organic layers and discard the aqueous layer. Wash the organic layer three times with 2-mL portions of saturated aqueous sodium chloride solution. Transfer the organic phase into a clean, dry Erlenmeyer flask and dry over granular anhydrous sodium sulfate for 15 minutes.

Sublimation Procedure. Set up a sublimation apparatus as shown in Technique 17, Figure 17.2C. (Alternatively, you can replace the sidearm test tube with a filter flask and the neoprene adapter can be replaced by a one-hole rubber stopper.) The sidearm test tube should be clean and dry. Transfer the dried organic layer to the sidearm test tube. In the hood, evaporate the solvent using a light stream of air while heating the tube in a warm (40°C) water bath. **Caution:** Excessive heating can lead to sublimation and product loss, so it is important to monitor the concentration step carefully, especially near the end when the solid product begins to form. You may find it useful to rotate the test tube in the air stream to remove the last traces of solvent. As soon as it is clear that there is no more solvent present, stop heating the sidearm test tube.

The technique for conducting a successful sublimation (Technique 17, Section 17.6B) is given below. *Your instructor should demonstrate the sublimation procedure before you start your sublimation.*

1. Make sure that all flammable solvents are removed from your workspace in the hood.
2. Assemble the apparatus shown in Technique 17, Figure 17.2C, using a sidearm test tube (or filter flask). Adjust the inner tube so that the bottom of the tube is about 1 cm above the solid in the sidearm test tube. Don't add water to the tube until Step 6.
3. Turn on the gas to a microburner and light it.
4. Make sure that the house vacuum or aspirator will pull a good vacuum. If the pressure is not sufficiently low, move to another place in the laboratory until you find a suitable place with a good vacuum. This is critical.
5. Attach the apparatus as shown in Figure 17.2C to the house vacuum or aspirator using heavy-walled pressure tubing, and turn the vacuum on.
6. Using a Pasteur pipette transfer ice water from an ice bath to the inner tube. The tube must be filled almost to the top, but don't allow water to flow over the top of the tube.
7. Hold the microburner by its base, and allow the flame to move around the bottom of the sidearm test tube. Do not hold the flame in one place! You should notice a white solid collecting on the inner tube. This is the purified camphor. Usually the material will sublime without melting, but if it does melt do not be concerned. You may need to heat the sides of the sidearm test tube to sublime all of the camphor. The sublimation process should not take more than 2 to 3 minutes. If it takes more than 3 minutes, you may need to remove the water in the inner tube and replace with new ice-cold water.

8. When all of the camphor has sublimed, turn off the microburner and allow the apparatus to cool. Do not turn off the vacuum! Remove the cold water from the central tube with a Pasteur pipette. When the apparatus has warmed up to room temperature, turn off the vacuum and disconnect the vacuum pressure tubing from the apparatus.
9. Carefully remove the inner tube and neoprene adapter (or one-hole rubber stopper) so that the purified camphor will not be dislodged from the inner tube.
10. Scrape off the camphor onto a piece of weighing paper. Camphor is a white solid. Determine the weight of the purified camphor and determine the percentage yield.
11. Determine the infrared spectrum of the sublimed camphor and compare to the spectrum provided in the experiment (use the dry film method in Technique 25, Section 25.4, or other method recommended by your instructor). The spectrum should demonstrate complete oxidation to camphor (absence of OH peak and presence of C=O peak). There should be sufficient material for the reduction of camphor to isoborneol, Part B. At the option of the instructor, determine the ^1H and ^{13}C NMR spectra of your camphor. Also, at the option of the instructor, determine the melting point (literature mp about 177°C , but it is often lower than this value).² Store the camphor in a tightly sealed vial.

Part B. Reduction of Camphor to Isoborneol.

Reductions. The camphor obtained in Part A should not contain borneol. If it does, show your infrared spectrum to your instructor and ask for advice. If the amount of camphor obtained in Part A is less than 0.1 g, obtain some camphor from the supply shelf to supplement your yield. If the amount is more than 0.1 g, scale up the reagents appropriately from the following amounts. Add 2.0 mL of methanol to the camphor contained in a 25-mL Erlenmeyer flask. Stir with a glass stirring rod until the camphor has dissolved. In portions, cautiously and intermittently add 0.10 g of sodium borohydride to the solution with a spatula. When all of the borohydride is added, boil the contents of the flask on a warm hot plate (low setting) for 2 minutes. Add more methanol if necessary to replace the solvent lost by evaporation.

Isolation and Analysis of Product. Allow the reaction mixture to cool for several minutes, and carefully add 5 mL of ice water. Collect the white solid by filtering it on a Hirsch funnel and, by using suction, allow the solid to dry for a few minutes. Transfer the solid to a dry Erlenmeyer flask. Add about 5 mL of methylene chloride to dissolve the product. Once the product has dissolved (add more solvent, if necessary), dry the solution over granular anhydrous sodium sulfate (see Technique 12, Section 12.9). When dry, the solution should not be cloudy. If the solution is still cloudy, add some more granular anhydrous sodium sulfate. Transfer the solution from the drying agent into a preweighed dry flask. Evaporate the solvent in a hood, as described in Part A.


²The observed melting point of camphor is often low. A small amount of impurity drastically reduces the melting point and increases the range (see Question 4 at the end of this experiment).

Determine the weight of the product, and calculate the percentage yield. If your instructor requests it, determine the melting point; pure isoborneol melts at 212°C. Determine the infrared spectrum of the product by the dry film method used previously with camphor. Compare the spectrum with the infrared spectra for borneol and isoborneol shown in the figures. Look especially to see if camphor has been completely reduced (absence of the C=O stretching band).


Part C. Percentages of Isoborneol and Borneol Obtained from the Reduction of Camphor.

NMR Determination. The percentage of each of the isomeric alcohols in the borohydride mixture can be determined from the NMR spectrum (see Technique 26). The 300-MHz NMR spectrum of pure borneol is shown near the end of this experiment. Notice that the hydrogen atom attached to the carbon bearing the hydroxyl group (—CH—OH) appears downfield at about 4.0 ppm due to the electronegative effect of the oxygen atom (see Technique 26, Table 26.1). In the 300-MHz NMR spectrum of pure isoborneol (near the end of this experiment), the corresponding hydrogen atom appears at 3.6 ppm. As shown on the last spectrum in this experiment, the 300-MHz proton-NMR spectrum of borohydride reduction product can be used to determine the percentages of borneol and isoborneol in a typical product obtained from reduction of camphor with sodium borohydride. Use the expansion of the region between 3.5 and 4.1 ppm and the individual integrals for the hydrogen atoms on the carbon bearing the OH groups. For this spectrum, the value for the integral rise is about 4.0 mm for borneol and 32.0 mm for isoborneol. As shown on the spectrum, we can use these integral values to calculate that this sample contains 89% isoborneol and 11% borneol. The calculated values for your product will likely vary from these values. Note that with many NMR spectrometers, the integration values of the signals will be printed next to or below the respective signals in the NMR spectrum. These values can then be used to determine the relative percentages of isoborneol and borneol in the product mixture.


Gas Chromatography. The isomer ratio and percentages can also be obtained by gas chromatography. Your instructor will provide instructions for preparing your sample. A Gow-Mac 69-360 instrument fitted with an 8-ft column of 10% Carbowax 20M, in an oven set at 180°C, and with a 40 mL/min helium flow rate will completely separate isoborneol and borneol from each other. In addition, any residual camphor can be observed. The retention times for camphor, isoborneol, and borneol are 8, 10, and 11 minutes, respectively. Other instrument conditions are provided in the Notes to the Instructor.

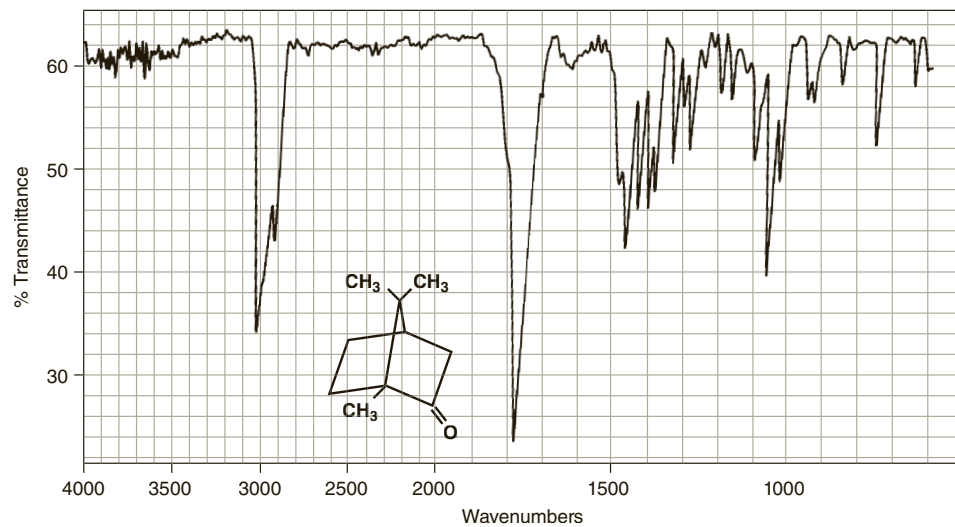


Infrared spectrum of camphor (KBr pellet).


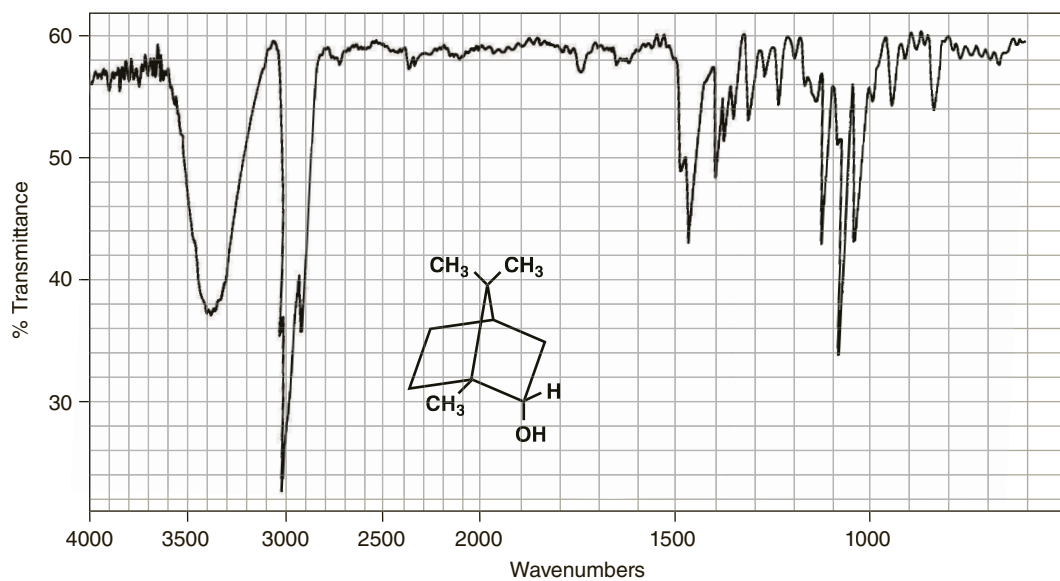


Infrared spectrum of borneol (KBr pellet).


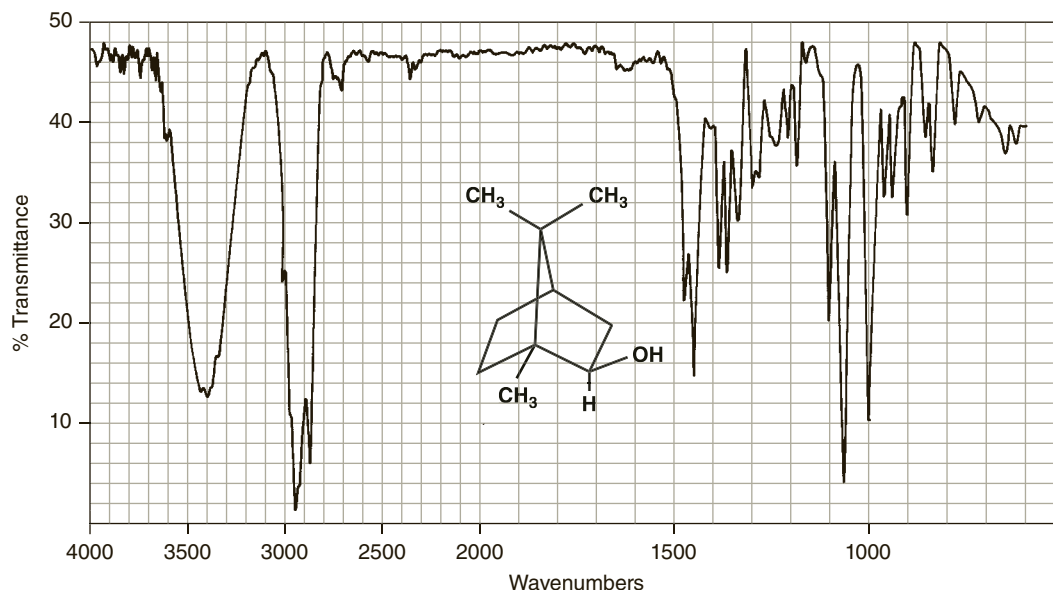


Infrared spectrum of isborneol (KBr pellet).


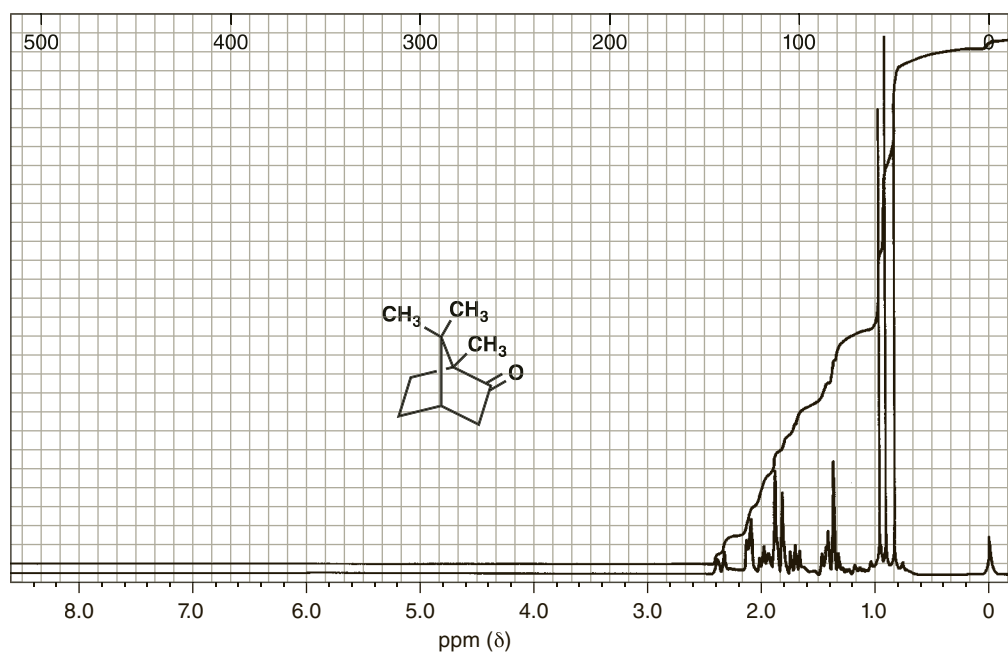
300-MHz NMR spectrum of camphor, CDCl₃.

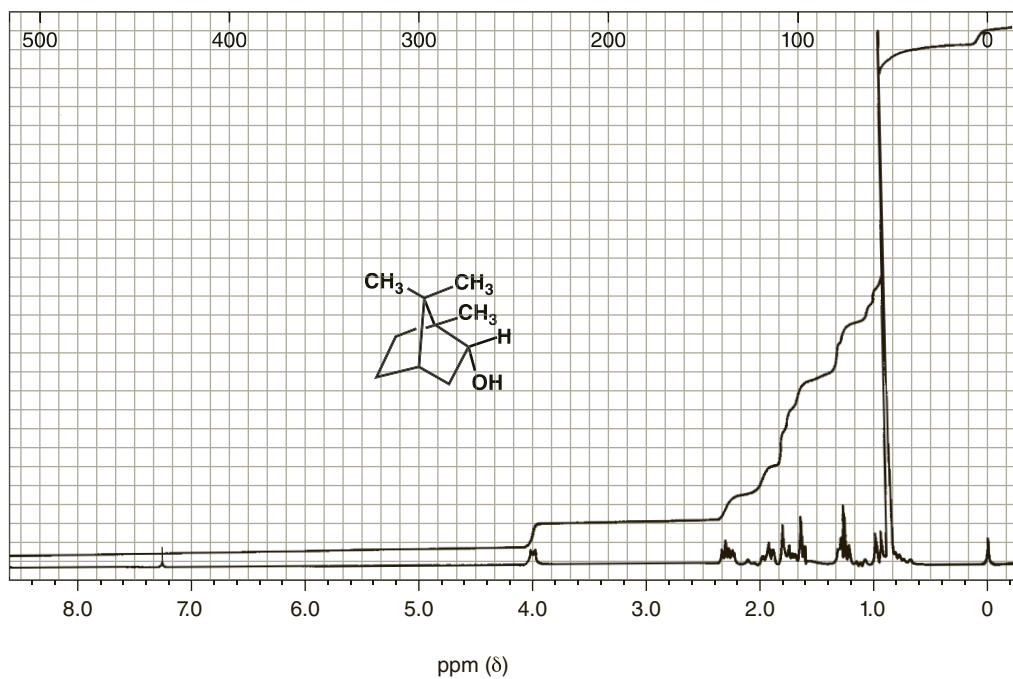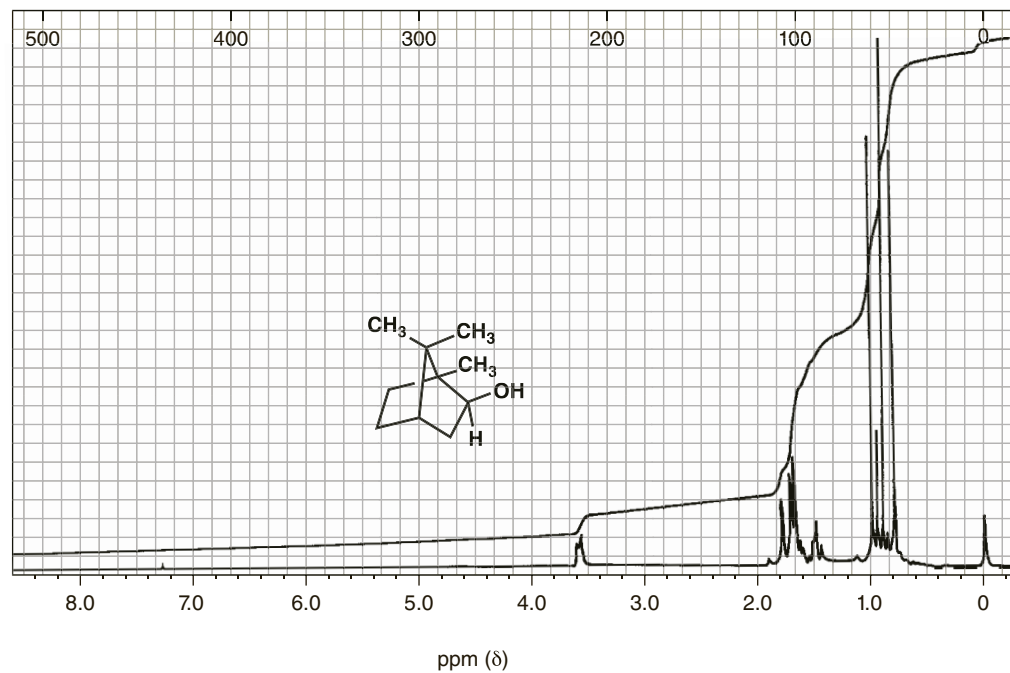300-MHz NMR spectrum of borneol, CDCl_3 .300-MHz NMR spectrum of isborneol, CDCl_3 .


a = 9.1 ppm q
b = 19.0 q
c = 19.6 q
d = 26.9 t
e = 29.8 t
f = 43.1 t
g = 43.1 d
h = 46.6 s
i = 57.4 s
j = 218.4 (not shown)

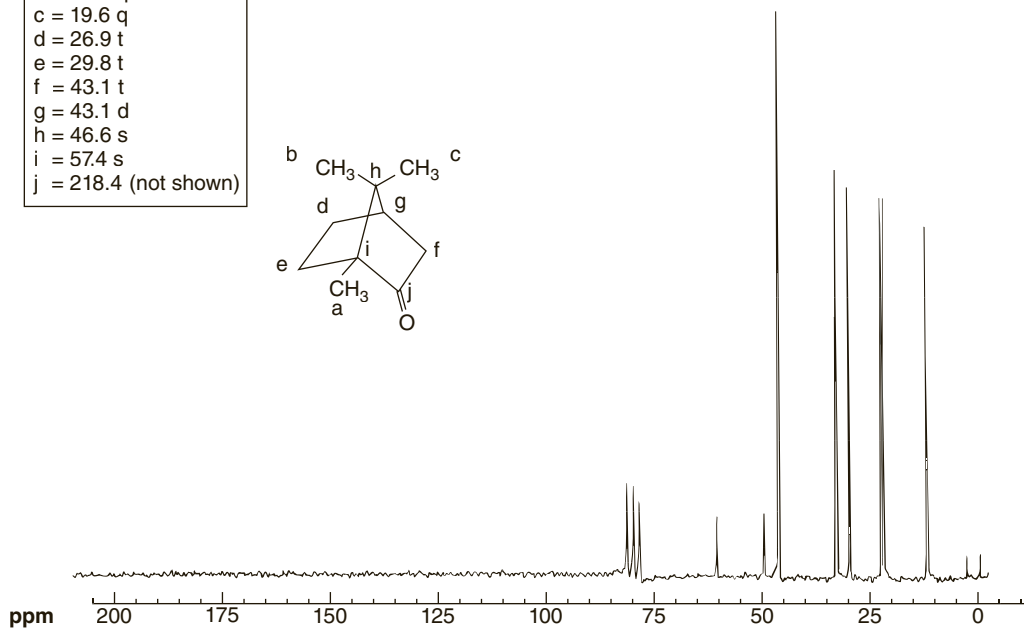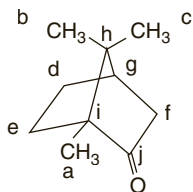

Carbon-13 spectrum of camphor, CDCl_3 .

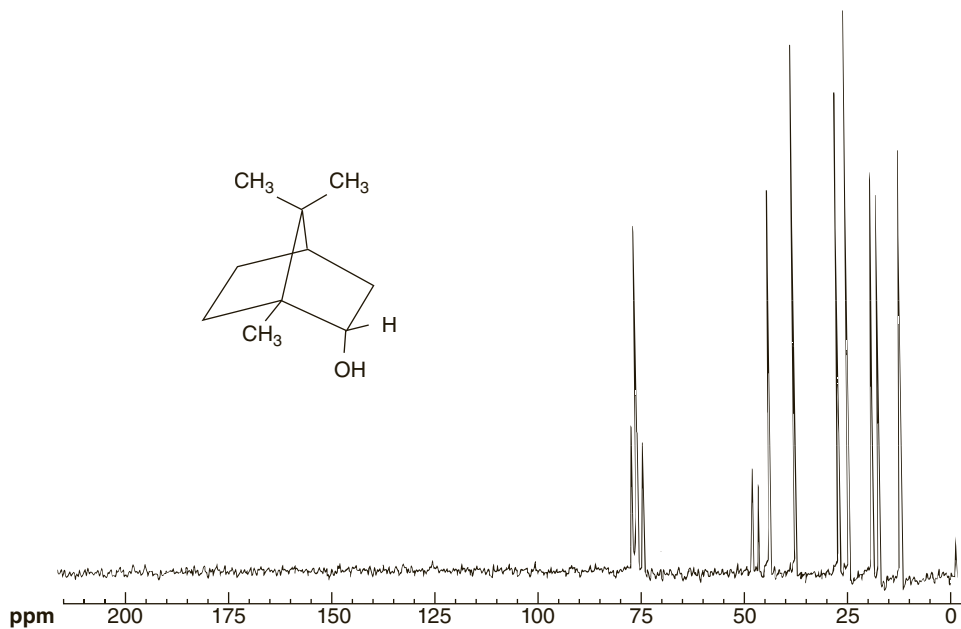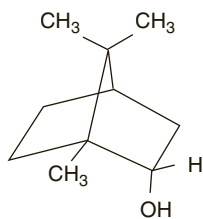

Carbon-13 spectrum of borneol, CDCl_3 .

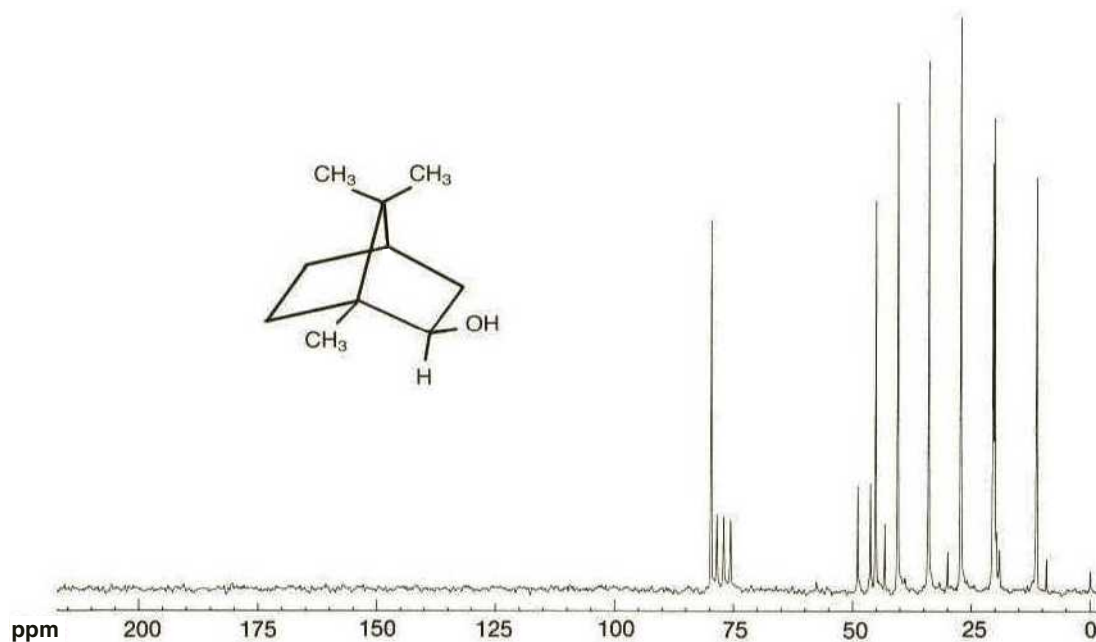


Carbon-13 spectrum of isborneol, CDCl_3 . (Small peaks at 9, 19, 30, and 43 are due to impurities.)

MOLECULAR MODELING (OPTIONAL)

In this exercise, we will seek to understand the experimental results obtained in the borohydride reduction of camphor and compare them to the results for the simpler norbornanone system (no methyl groups). Because the hydride ion is an electron donor, it must place its electrons into an empty substrate orbital to form a new bond. The most logical orbital for this action is the LUMO (lowest unoccupied molecular orbital). Accordingly, the focus of our calculations will be the shape and location of the LUMO.

Part A.

Build a model of norbornanone, and submit it to an AM1-level calculation of its energy using a geometry optimization. Also request that density and LUMO surfaces be calculated, along with a density–LUMO surface (a mapping of the LUMO onto the density surface).

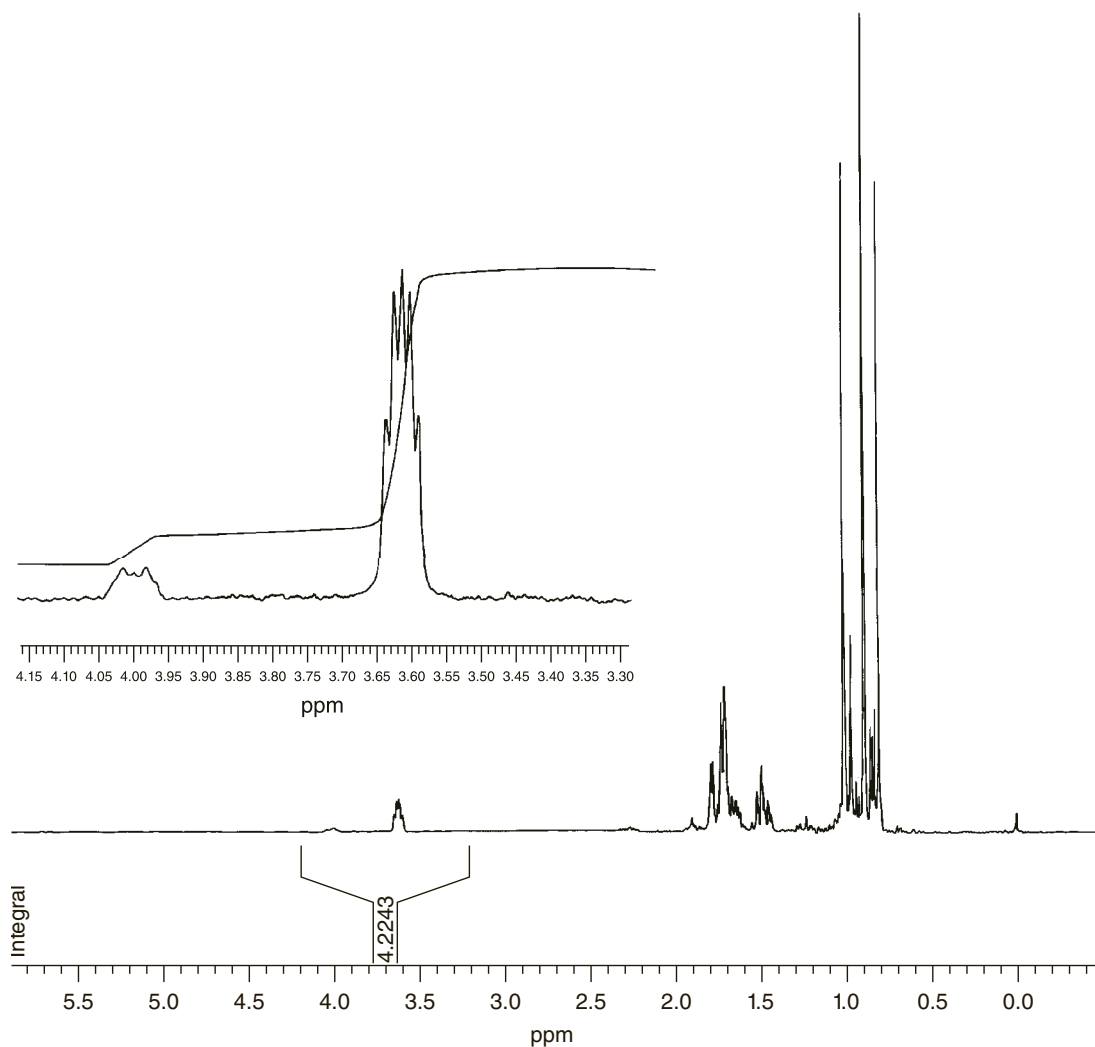

300-MHz Proton NMR spectrum of borohydride reduction product, CDCl_3 . Inset: Expansion of the 3.5–4.1 ppm region.

When the calculation is complete, display the LUMO on the norbornanone skeleton. Where is the size of the LUMO (its density) largest? Which atom is this? This is the expected site of addition. Now map a density surface onto the same norbornanone surface. When you consider the approach of the borohydride ion, which face is less hindered? Is an *exo* or *endo* approach favored? An easier way to decide is to view the density–LUMO surface. On this surface, the intersection of the LUMO with the density surface is color-coded. The spot where the access to the LUMO is easiest (the location of its largest value) will be coded blue. Is this spot on the *endo* or on the *exo* face? Do your modeling results agree with the observed reaction percentages (see Part C completed earlier)?

Part B.

Follow the same instructions given earlier for norbornanone using camphor—that is, calculate and view density, LUMO, and density–LUMO surfaces.

Do you reach the same conclusions as for norbornanone? Are there new stereochemical considerations? Do your conclusions agree with the experimental results (the borneol/isborneol ratio) you obtained in this experiment? In your report, discuss your modeling results and how they relate to your experimental results.

REFERENCES

- Brown, H. C.; Muzzio, J. Rates of Reaction of Sodium Borohydride with Bicyclic Ketones. *J. Am. Chem. Soc.* **1966**, *88*, 2811.
- Dauben, W. G.; Fonken, G. J.; Noyce, D. S. Stereochemistry of Hydride Reductions. *J. Am. Chem. Soc.* **1956**, *78*, 2579.
- Lang, P. T.; Harned, A. M.; Wissinger, J. E. Oxidation of Borneol to Camphor Using Oxone® and Catalytic Sodium Chloride: A Green Experiment for the Undergraduate of Organic Chemistry Laboratory. *J. Chem. Educ.* **2011**, *88*, 652.
- Markgraf, J. H. Stereochemical Correlations in the Camphor Series. *J. Chem. Educ.* **1967**, *44*, 36.

QUESTIONS

1. Interpret the major absorption bands in the infrared spectra of camphor, borneol, and isoborneol.
2. Explain why the *gem*-dimethyl groups appear as separate peaks in the proton-NMR spectrum of isoborneol, although they almost overlap in borneol.
3. A sample of isoborneol prepared by reduction of camphor was analyzed by infrared spectroscopy and showed a band at 1750 cm^{-1} . This result was unexpected. Why?
4. The observed melting point of camphor is often low. Look up the molal freezing point–depression constant *K* for camphor, and calculate the expected depression of the melting point of a quantity of camphor that contains 0.5 molal impurity. (*Hint*: Look in a general chemistry book under “freezing-point depression” or “colligative properties of solutions.”)
5. The peak assignments are shown on the carbon-13 NMR spectrum of camphor. Using these assignments as a guide, assign as many peaks as possible in the carbon-13 spectra of borneol and isoborneol.

Multistep Reaction Sequences: The Conversion of Benzaldehyde to Benzilic Acid

Green chemistry

Multistep reactions

Thiamine-catalyzed reaction

Oxidation with nitric acid

Rearrangement

Crystallization

Computational chemistry (optional)

The experiment demonstrates multistep synthesis of benzilic acid starting from benzaldehyde. In Experiment 34A, benzaldehyde is converted to benzoin using a thiamine-catalyzed reaction. This part of the experiment demonstrates how a “green” reagent can be utilized in organic chemistry. In Experiment 34B, nitric acid oxidizes benzoin to benzil. Finally, in Experiment 34C, benzil is rearranged to benzilic acid. The scheme on the next page shows the reactions.


REQUIRED READING

<i>Review:</i> Technique 6	Heating and Cooling Methods, Sections 6.1–6.3
Technique 7	Reaction Methods, Sections 7.1–7.4
Technique 8	Filtration, Section 8.3
Technique 9	Physical Constants of Solids: The Melting Point, Sections 9.7 and 9.8
Technique 11	Crystallization: Purification of Solids, Section 11.3
Technique 12	Extractions, Separations, and Drying Agents, Section 12.4
Technique 25	Infrared Spectroscopy, Section 25.4
<i>New:</i> Essay and Experiment 20	Computational Chemistry (Optional)

NOTES TO THE INSTRUCTOR

Although this experiment is intended to illustrate a multistep synthesis to the students, each part may be done separately, or two out of the three reactions can


be linked together. The sections on “Special Instructions” and “Suggested Waste Disposal” are included in each part of this experiment.

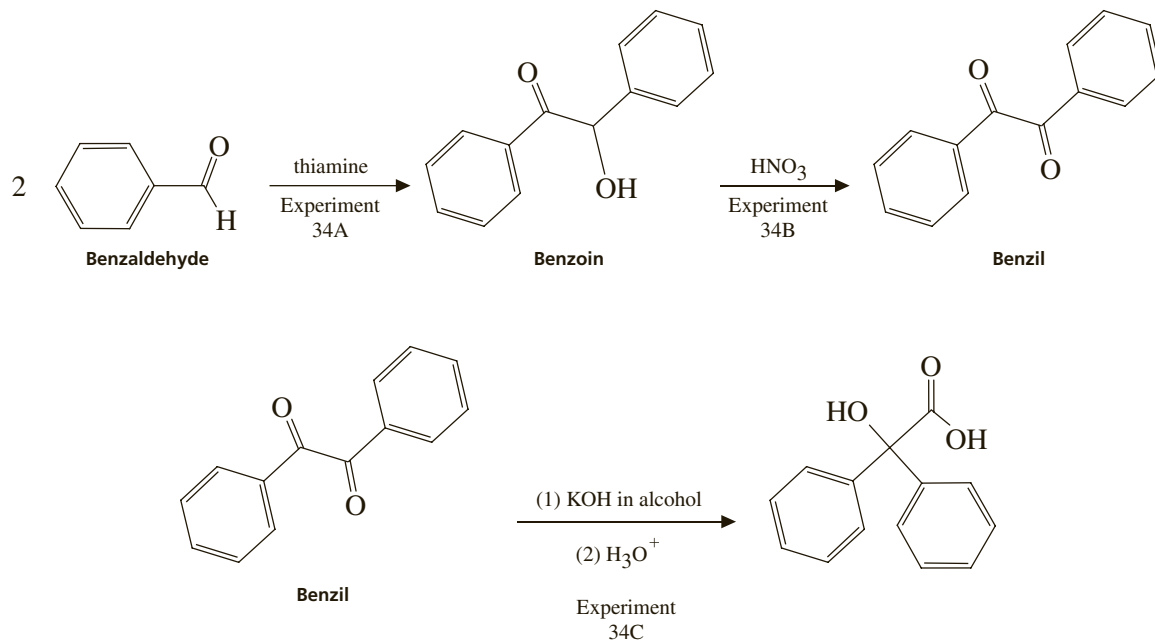


34A EXPERIMENT 34A

Preparation of Benzoin by Thiamine Catalysis

In this experiment, two molecules of benzaldehyde will be converted to benzoin using the catalyst thiamine hydrochloride. This reaction is known as a benzoin condensation reaction:

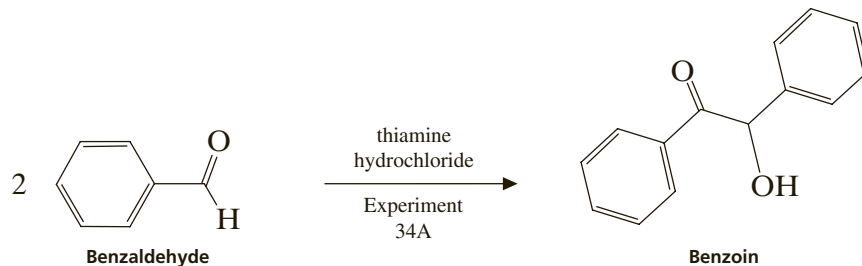


Thiamine hydrochloride is structurally similar to thiamine pyrophosphate (TPP). TPP is a coenzyme universally present in all living systems. It catalyzes several biochemical reactions in natural systems. It was originally discovered as a required nutritional factor (vitamin) in humans by its link with the disease beriberi. **Beriberi** is a disease of the peripheral nervous system caused by a deficiency of Vitamin B₁ in the diet. Symptoms include pain and paralysis of the extremities, emaciation, and swelling of the body. The disease is most common in Asia.

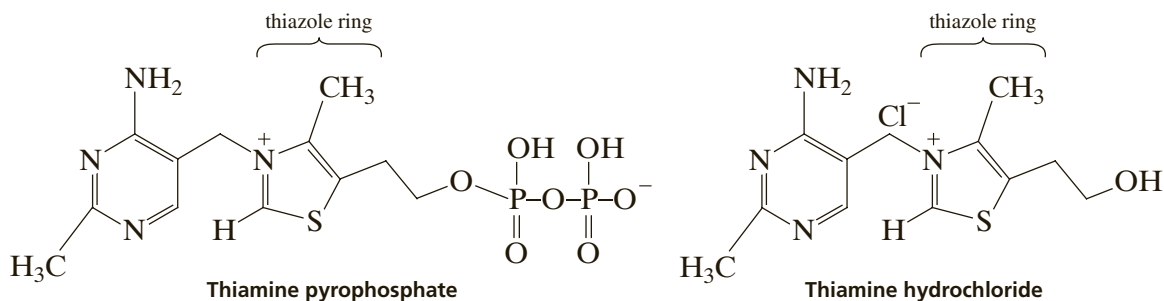


Thiamine binds to an enzyme before the enzyme is activated. The enzyme also binds to the substrate (a large protein). Without the coenzyme thiamine, no chemical reaction would occur. The coenzyme is the **chemical reagent**. The protein molecule (the enzyme) helps and mediates the reaction by controlling stereochemical, energetic, and entropic factors, but it is nonessential to the overall course of reactions that it catalyzes. A special name, vitamins, is given to coenzymes that are essential to the nutrition of the organism.

The most important part of the entire thiamine molecule is the central ring, the thiazole ring, which contains nitrogen and sulfur. This ring constitutes the **reagent** portion of the coenzyme. Experiments with the model compound 3,4-dimethylthiazolium bromide have explained how thiamine-catalyzed reactions work. It was found that this model thiazolium compound rapidly exchanged the C-2 proton for deuterium in D_2O solution. At a pD of 7 (no pH here), this proton was completely exchanged in seconds!


This indicates that the C-2 proton is more acidic than one would have expected. It is apparently easily removed because the conjugate base is a highly stabilized ylide. An **ylide** is a compound or intermediate with positive and negative formal charges on adjacent atoms.

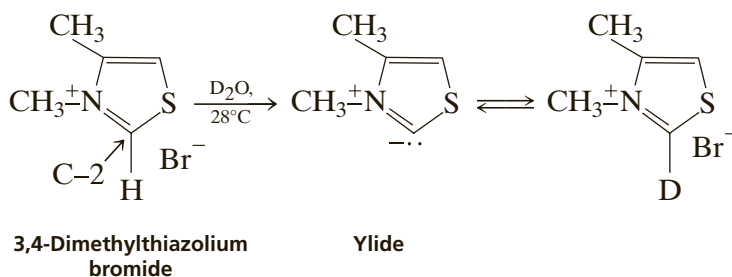


The sulfur atom plays an important role in stabilizing this ylide. This was shown by comparing the rate of exchange of 1,3-dimethylimidazolium ion with the rate for the thiazolium compound shown in the previous equation. The dinitrogen compound exchanged its C-2 proton more slowly than the sulfur-containing ion. Sulfur, being in the third row of the periodic chart, has *d* orbitals available for bonding to adjacent atoms. Thus, it has fewer geometrical restrictions than carbon and nitrogen atoms do and can form carbon–sulfur multiple bonds in situations in which carbon and nitrogen normally would not.

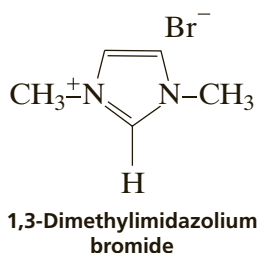

In Experiment 34A, we will utilize thiamine hydrochloride rather than TPP to catalyze the benzoin condensation. The mechanism is shown on the next page. For simplicity, only the thiazole ring is shown.

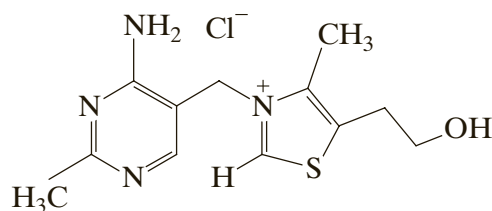

Thiamine hydrochloride

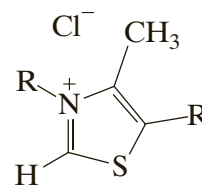

Thiazole ring in thiamine hydrochloride


The mechanism involves the removal of the proton at C-2 from the thiazole ring with a weak base to give the ylide (Step 1). The ylide acts as a nucleophile that adds to the carbonyl group of benzaldehyde forming an intermediate (Step 2). A proton is removed to yield a new intermediate with a double bond (Step 3). Notice that the nitrogen atom helps to increase the acidity of that proton. This intermediate can now react with a second benzaldehyde molecule to yield a new intermediate (Step 4). A base removes a proton to produce benzoin and also regenerates the ylide (Step 5). The ylide reenters the mechanism to form more benzoin by the condensation of two more molecules of benzaldehyde.

SPECIAL INSTRUCTIONS

This experiment may be conducted concurrently with another experiment. It involves a few minutes at the beginning of a laboratory period for mixing reagents. The remaining portion of the period may be used for another experiment.

SUGGESTED WASTE DISPOSAL

Pour all of the aqueous solutions produced in this experiment into a waste container designated for aqueous waste. The ethanolic mixtures obtained from the crystallization of crude benzoin should be poured into a waste container designated for nonhalogenated waste.

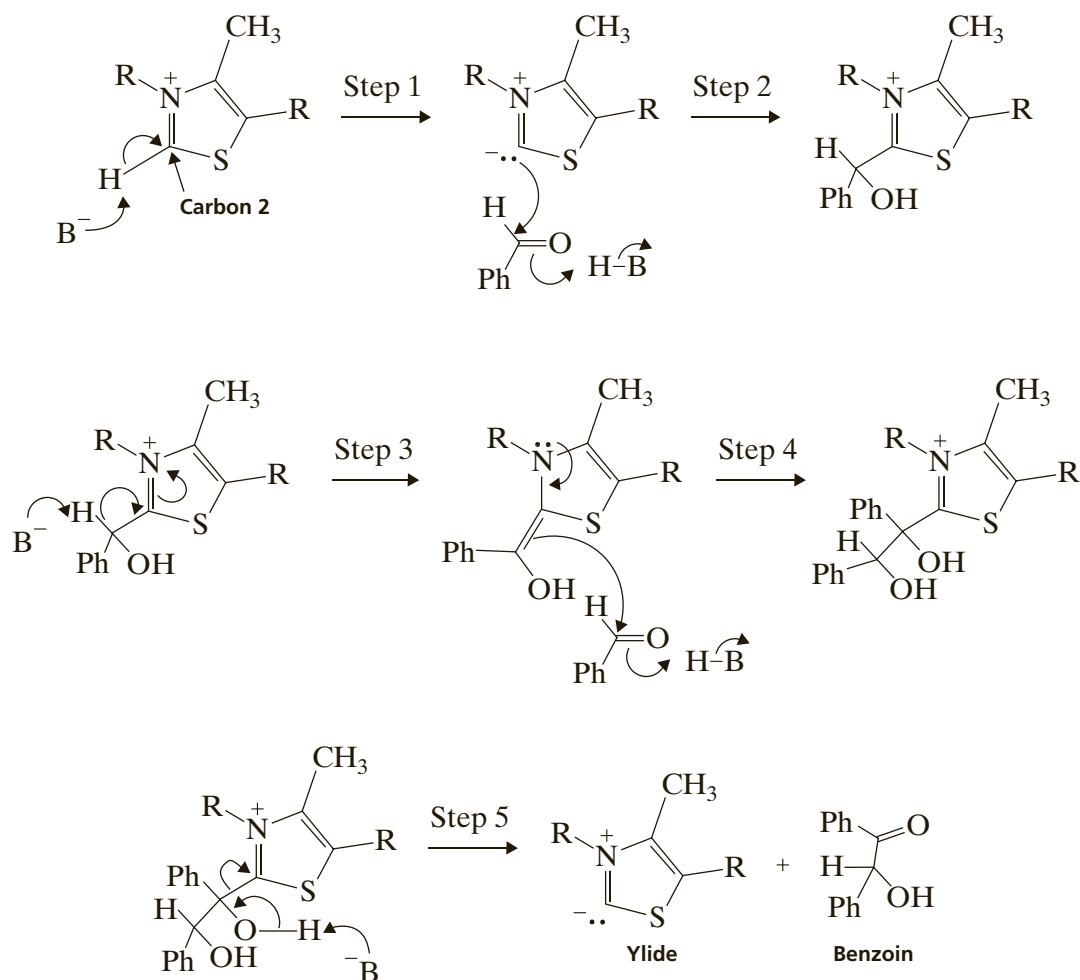

NOTES TO THE INSTRUCTOR

It is essential that the benzaldehyde used in this experiment be *pure*. Benzaldehyde is easily oxidized in air to benzoic acid. Even when benzaldehyde *appears* free of benzoic acid by infrared spectroscopy, you should check the purity of your benzaldehyde and thiamine by following the instructions given in the first paragraph of the Procedure (“Reaction Mixture”). When the benzaldehyde is pure, the solution will be nearly filled with solid benzoin after 2 days (you may need to scratch the inside of the flask to induce crystallization). If no solid appears or very little appears, then there is a problem with the purity of the benzaldehyde. If possible, use a newly opened bottle that has been purchased recently. However, it is essential that you check both the old and new benzaldehyde before doing the laboratory experiment.

We have found that the following procedure does an adequate job of purifying benzaldehyde. The procedure does not require distillation of benzaldehyde. Shake the benzaldehyde in a separatory funnel with an equal volume of 5% aqueous sodium carbonate solution. Shake gently and occasionally open the stopcock of the funnel to vent carbon dioxide gas. An emulsion forms that may take 2–3 hours to separate. It is helpful to stir the mixture occasionally during this period to help break the emulsion. Remove the lower sodium carbonate layer, including any remaining emulsion. Add about $\frac{1}{4}$ volume of water to the benzaldehyde and shake the mixture gently to avoid an emulsion. Remove the cloudy lower organic layer and dry the benzaldehyde with calcium chloride until the next day. Any remaining

cloudiness is removed by gravity filtration through fluted filter paper. The resulting clear, purified benzaldehyde should be suitable for this experiment without vacuum distillation. You *must* check the purified benzaldehyde to see if it is suitable for the experiment by following the instructions in the first paragraph of the Procedure.

It is advisable to use a fresh bottle of thiamine hydrochloride, which should be stored in the refrigerator. Fresh thiamine does not seem to be as important as pure benzaldehyde for success in this experiment.

PROCEDURE

Reaction Mixture

Add 0.30 g thiamine hydrochloride to a 25-mL Erlenmeyer flask. Dissolve the solid in 0.45 mL of water by swirling the flask. Add 3.0 mL of 95% ethanol and swirl the solution until it is homogeneous. To this solution, add 0.90 mL of an aqueous sodium hydroxide¹ solution and swirl the flask until the bright yellow color fades to a pale yellow. Weigh the flask and solution, add 0.90 mL of benzaldehyde, and reweigh the flask to determine an accurate weight of benzaldehyde introduced to the flask. Swirl the contents of the flask until it is homogeneous. Stopper the flask and let it stand in a dark place for at least 2 days.

Isolation of Crude Benzoin

If crystals have not formed after 2 days, initiate crystallization by scratching the inside of the flask with a glass stirring rod. Allow about 5 minutes for the crystals of benzoin to form fully. Place the flask, with crystals, into an ice bath for 5–10 minutes.

If for some reason the product separates as an oil, it may be helpful to scratch the flask with a glass rod or seed the mixture by allowing a small amount of solution to dry on the end of a glass rod and then placing this into the mixture. Cool the mixture in an ice bath before filtering.

Break up the crystalline mass with a spatula, swirl the flask rapidly, and quickly transfer the benzoin to a Hirsch funnel under vacuum (Technique 8, Section 8.3, and Figure 8.5). Wash the crystals with three 1-mL portions of ice-cold water. Allow the benzoin to dry in the Hirsch funnel by drawing air through the crystals for about 5 minutes. Transfer the benzoin to a watch glass and allow it to dry in air until the next laboratory period. The product may also be dried in a few minutes in an oven set at about 100°C.

Yield Calculation and Melting-Point Determination

Weigh the benzoin and calculate the percentage yield based on the amount of benzaldehyde used initially. Determine the melting point (pure benzoin melts between 134°C and 135°C). Because your crude benzoin will normally melt between 129°C and 132°C, the benzoin should be crystallized before the conversion to benzil (Experiment 34B).

Crystallization of Benzoin

Purify the crude benzoin by crystallization from hot 95% ethanol (use 0.8 mL of alcohol/0.1 g of crude benzoin) using a 10-mL Erlenmeyer flask for the crystallization (Technique 11, Section 11.3, omit step 2 shown in Figure 11.4). After the crystals

¹ Dissolve 8.0 g of NaOH in 100 mL water.


have cooled in an ice bath, collect the crystals on a Hirsch funnel. The product may be dried in a few minutes in an oven set at about 100°C. Determine the melting point of the purified benzoin. If you are not scheduled to perform Experiment 34B, submit the sample of benzoin, along with your report, to the instructor.

Spectroscopy


Determine the infrared spectrum of the benzoin by the dry film method (Technique 25, Section 25.4). A spectrum is shown here for comparison.

QUESTIONS

1. The infrared spectra of benzoin and benzaldehyde are given in this experiment. Interpret the principal peaks in the spectra.

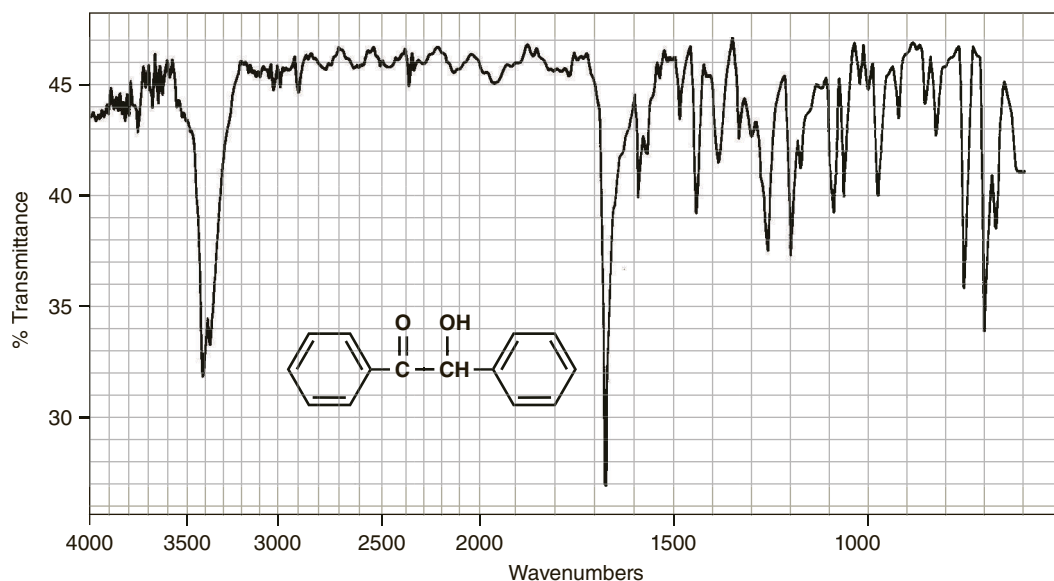


Infrared spectrum of benzoin, KBr.

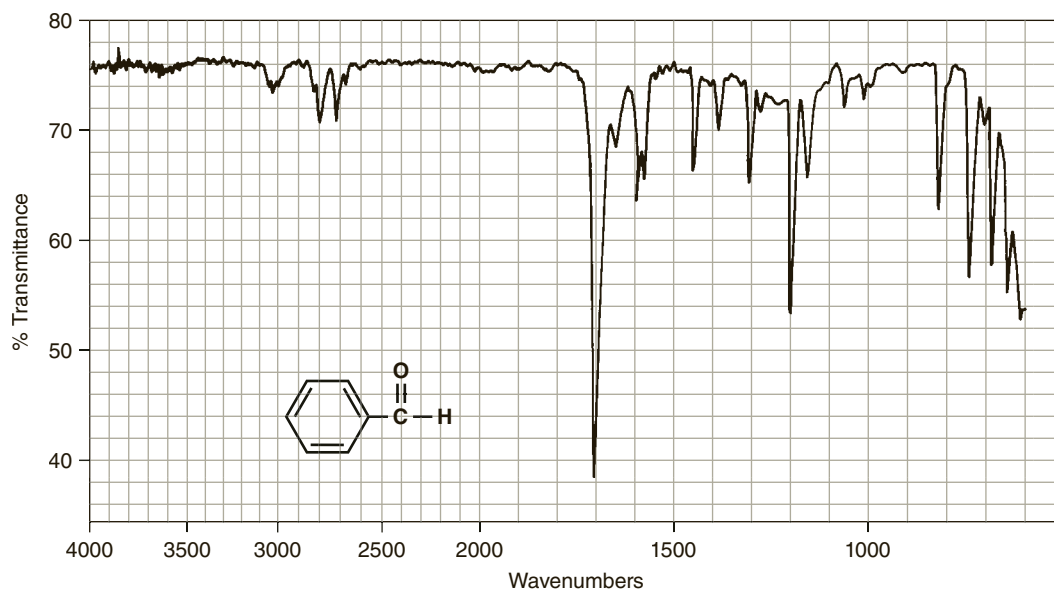

Infrared spectrum of benzaldehyde (neat).


- How do you think the appropriate enzyme would have affected the reaction (degree of completion, yield, stereochemistry)?
- What modifications of conditions would be appropriate if the enzyme were to be used?
- Draw a mechanism for the cyanide-catalyzed conversion of benzaldehyde to benzoin. The intermediate, shown in brackets, is thought to be involved in the mechanism.

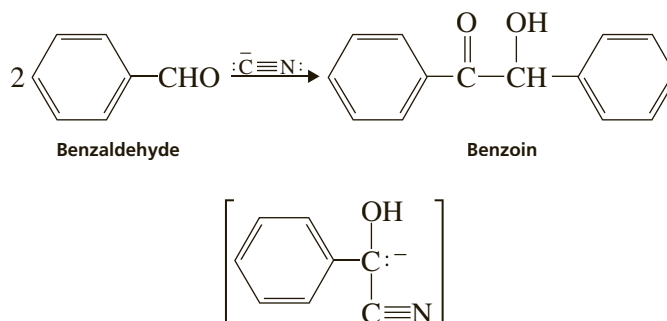

34B EXPERIMENT 34B

Preparation of Benzil

In this experiment, benzil is prepared by the oxidation of an α -hydroxyketone, benzoin. This experiment uses the benzoin prepared in Experiment 34A and is the second step in the multistep synthesis. This oxidation can be done easily with mild oxidizing agents such as Fehling's solution (alkaline cupric tartrate complex) or copper sulfate in pyridine. In this experiment, the oxidation is performed with nitric acid.

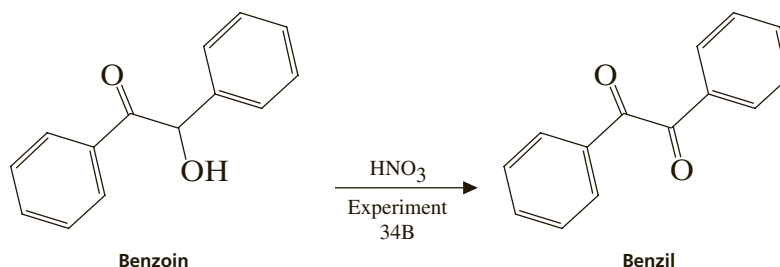

SPECIAL INSTRUCTIONS

Nitric acid should be dispensed in a good hood to avoid the choking odor of this substance. The vapors will irritate your eyes. Avoid contact with your skin. During the reaction, considerable amounts of noxious nitrogen oxide gases are evolved. Be sure to run the reaction in a good fume hood.

SUGGESTED WASTE DISPOSAL

The aqueous nitric acid wastes should be poured into a container designated for nitric acid wastes. Do not put them into the aqueous waste container. The ethanolic wastes from the crystallization should be poured into the nonhalogenated waste container.

PROCEDURE

Reaction Mixture


Place 0.30 g of benzoin (Experiment 34A) into a 5-mL conical vial and add 1.5 mL of concentrated nitric acid. Add a magnetic spin vane and attach an air condenser. In a hood, set up the apparatus for heating in a hot water bath, as shown in Technique 6, Figure 6.6. Heat the mixture in a hot water bath at about 70°C for 1 hour, with stirring. Avoid heating the mixture above this temperature to reduce the possibility of forming a by-product.¹ During the 1-hour heating period, nitrogen oxide gases (red) will be evolved. If it appears that gases are still being evolved after 1 hour, continue heating for another 15 minutes but then discontinue heating at that time.

Isolation of Crude Benzil

Cool the mixture for a few minutes and detach the air condenser. With a Pasteur pipette, transfer the reaction mixture to a beaker containing 4 mL of ice-cold water. Rinse the conical vial and spin vane with a small amount of water. Cool the mixture in an ice bath until the crystals have formed. If the material oils out rather than crystallizes, scratch the oil vigorously with a spatula until it does crystallize completely. Collect the crude product on a Hirsch funnel under vacuum. Wash it well with cold water (about 5 mL) to remove the nitric acid. Continue drawing air through the solid mass on the Hirsch funnel to help dry the solid. Weigh the solid.

Crystallization of Product

Purify the solid by dissolving it in hot 95% ethanol in a small Erlenmeyer flask (about 5 mL per 0.5 g of product) using a hot plate as the heating source. Be careful not to melt the solid on the hot plate. You can avoid melting the benzil by occasionally lifting the flask from the hot plate and swirling the contents of the flask. You want the solid to dissolve in the hot solvent rather than melt. You will obtain better

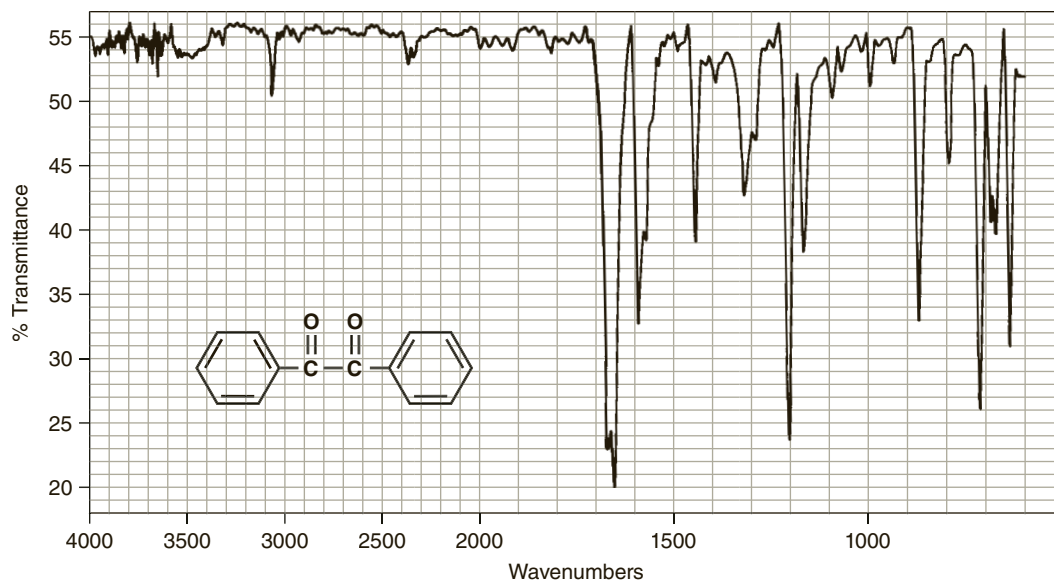

Infrared spectrum of benzil, KBr.

¹ At higher temperatures, some 4-nitrobenzil will be formed along with benzil.

crystals if you add a little extra solvent after it dissolves completely. Remove the flask from the hot plate and allow the solution to cool slowly. As the solution cools, seed it with a solid product that forms on a spatula after the spatula is dipped into the solution. The solution may become supersaturated unless this is done, and crystallization will occur too rapidly. Yellow crystals are formed. Cool the mixture in an ice bath to complete the crystallization. Collect the product on a Hirsch funnel under vacuum. Rinse the flask with small amounts (about 1 mL total) of ice-cold 95% ethanol to complete the transfer of product to the Hirsch funnel. Continue drawing air through the crystals on the Hirsch funnel by suction for about 5 minutes. Then remove the crystals and air-dry them.


Yield Calculation and Melting-Point Determination

Weigh the dry benzil and calculate the percentage yield. Determine the melting point. The melting point of pure benzil is 95°C. Submit the benzil to the instructor unless it is to be used to prepare benzilic acid (Experiment 34C). Obtain the infrared spectrum of benzil using the dry film method. Compare the spectrum to the one shown. Also compare the spectrum with that of benzoin. What differences do you notice?

34C EXPERIMENT 34C

Preparation of Benzilic Acid

In this experiment, benzilic acid will be prepared by causing the rearrangement of the α -diketone benzil. Preparation of benzil is described in Experiment 34B. The rearrangement of benzil proceeds in the following way:

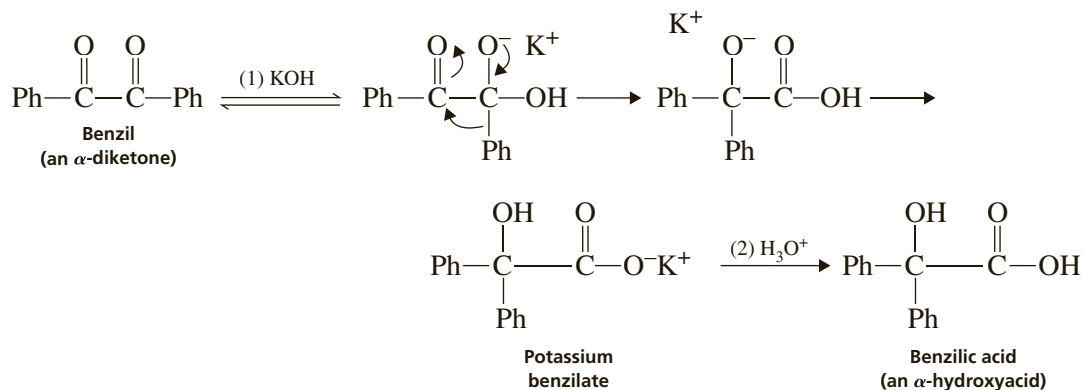

The driving force for the reaction is provided by the formation of a stable carboxylate salt (potassium benzilate). Once this salt is produced, acidification yields benzilic acid. The reaction can generally be used to convert aromatic α -diketones to aromatic α -hydroxyacids. Other compounds, however, also will undergo benzilic acid-type of rearrangement (see questions).

SPECIAL INSTRUCTIONS

This experiment works best with pure benzil. The benzil prepared in Experiment 34B is usually of sufficient purity after it has been crystallized.

SUGGESTED WASTE DISPOSAL

Pour all of the aqueous filtrates into the waste bottle designated for aqueous waste. Ethanolic filtrates should be put in the nonhalogenated organic waste bottle.

PROCEDURE

Running the Reaction

Add 0.100 g benzil and 0.30 mL 95% ethanol to a 3-mL conical vial. Place a spin vane in the vial and attach an air condenser. Heat the mixture with an aluminum block (90–100°C) while stirring until the benzil has dissolved (see inset in Technique 6, Figure 6.1A). Using a 9-inch Pasteur pipette, add dropwise 0.25 mL of an aqueous potassium hydroxide solution¹ downward through the condenser into the vial. Gently boil the mixture (aluminum block about 110°C) while stirring for 15 minutes. The mixture will be blue-black. As the reaction proceeds, the color will turn to brown, and the solid should dissolve completely. Solid potassium benzilate may form during the reaction period. At the end of the heating period, remove the assembly from the aluminum block and allow it to cool for a few minutes.

Crystallization of Potassium Benzilate

Detach the air condenser when the apparatus is cool enough to handle. Transfer the reaction mixture, which may contain some solid, with a Pasteur pipette into a 10-mL beaker. Allow the mixture to cool to room temperature and then cool in an ice-water bath for about 15 minutes until crystallization is complete. It may be necessary to scratch the inside of the beaker with a glass stirring rod to induce crystallization. Crystallization is complete when virtually the entire mixture has solidified. Collect the crystals on a Hirsch funnel by vacuum filtration and wash the crystals thoroughly with three 1-mL portions of ice-cold 95% ethanol. The solvent should remove most of the color from the crystals.

Transfer the solid, which is mainly potassium benzilate, to a 10-mL Erlenmeyer flask containing 3 mL of hot (70°C) water. Stir the mixture until all solid has dissolved or until it appears that the remaining solid will not dissolve. Any remaining solid will likely form a fine suspension. If solid does remain in the flask, filter the mixture in the following manner. Place about 0.5 g of Celite (Filter Aid) in a beaker with about 5 mL of water. Stir the mixture vigorously and then pour the contents into a Hirsch funnel (with filter paper) or a small Büchner funnel while applying a *gentle* vacuum, as in a vacuum filtration (Technique 8, Section 8.3, and Figure 8.5). Be careful not to let the Celite dry completely. This procedure will cause a thin layer of Celite to be deposited on the filter paper. Discard the water that passes through this filter. Pass the mixture containing potassium benzilate through this filter, using *gentle* suction. The filtrate should be clear. Transfer the filtrate to a 10-mL Erlenmeyer flask. If no solid remains in the flask, the filtration step may be omitted. In either case, proceed to the next step.


¹ The aqueous potassium hydroxide solution should be prepared for the class by dissolving 2.75 g of potassium hydroxide in 6.0 mL of water. This will provide enough solution for 20 students, assuming little solution is wasted.

Formation of Benzilic Acid

With swirling of the flask, add dropwise 0.50 mL of 1M hydrochloric acid to the solution of potassium benzilate. As the solution becomes acidic, solid benzilic acid will begin to precipitate. The pH should be about 2; if it is higher than this, add a few more drops of acid and check the pH again. Allow the mixture to cool to room temperature and then complete the cooling in an ice bath. Collect the benzilic acid by vacuum filtration using a Hirsch funnel. Wash the crystals with 3–4 mL of ice-cold water to remove potassium chloride salt that sometimes coprecipitates with benzilic acid during the neutralization with hydrochloric acid. Remove the wash water by drawing air through the filter. Dry the product thoroughly by allowing it to stand until the next laboratory period.

Melting Point and Crystallization of Benzilic Acid

Weigh the dry benzilic acid and determine the percentage yield. Determine the melting point of the dry product. Pure benzilic acid melts at 150°C. If necessary, crystallize the product from hot water using a Craig tube (Technique 11, Section 11.4, and Figure 11.6). If some impurities remain undissolved, filter the mixture using the following procedure. It will be necessary to keep the mixture hot during this filtration step. Transfer the hot mixture to a test tube with a Pasteur pipette. Clean the Craig tube and filter the mixture by transferring it back to the Craig tube with a filter-tip pipette. Cool the solution and induce crystallization, if necessary. Allow the mixture to stand at room temperature until crystallization is complete (about 15 minutes). Cool the mixture in an ice bath and collect the crystals by centrifugation. Determine the melting point of the crystallized product after it is thoroughly dry.

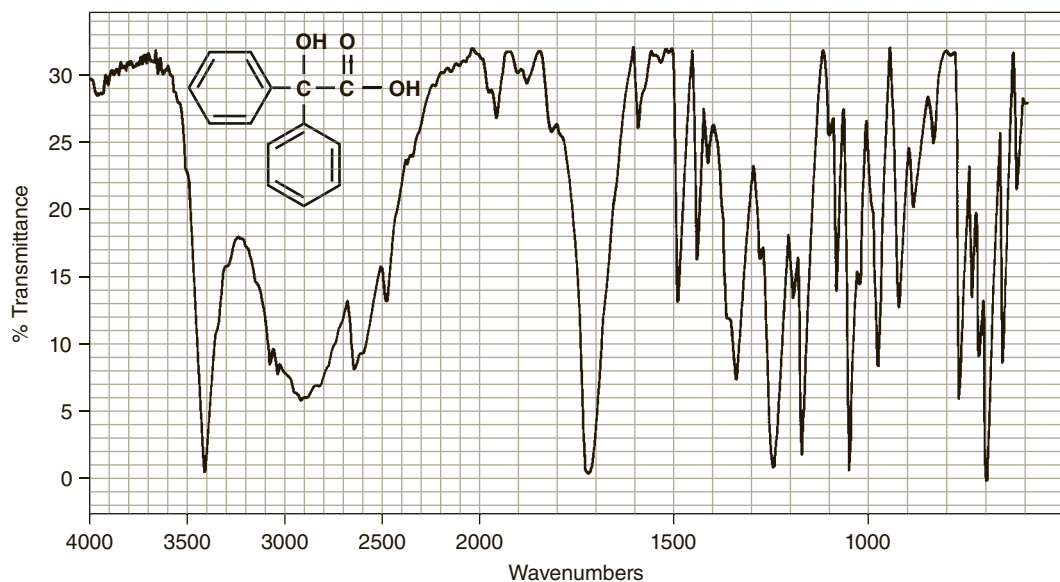


Infrared spectrum of benzilic acid, KBr.




At the instructor's option, determine the infrared spectrum of the benzilic acid in potassium bromide (Technique 25, Section 25.5). Calculate the percentage yield. Submit the sample to your laboratory instructor in a labeled vial.

QUESTIONS

1. Show how to prepare the following compounds, starting from the appropriate aldehyde.

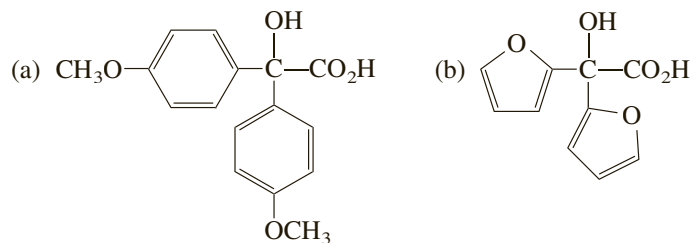

2. Give the mechanisms for the following transformations:

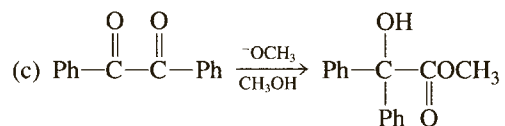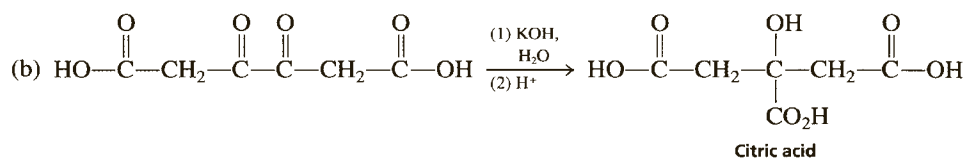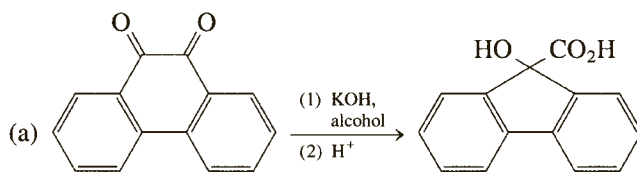

3. Interpret the infrared spectrum of benzoic acid.

Triphenylmethanol and Benzoic Acid


Grignard reactions

Extractions

Crystallization



Use of a separatory funnel


In this experiment, you will prepare a Grignard reagent, or organomagnesium reagent. The reagent is phenylmagnesium bromide.

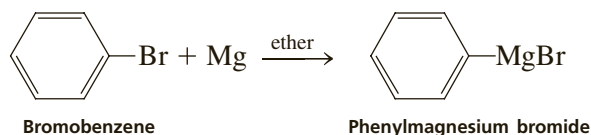


This reagent will be converted to a tertiary alcohol or a carboxylic acid, depending on the experiment selected.

Triphenylmethanol (Experiment 35A)

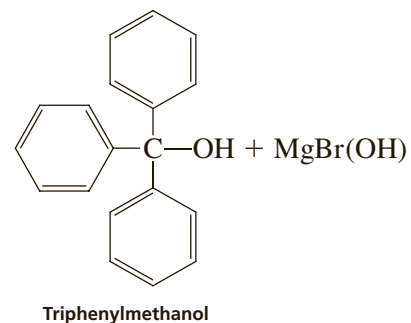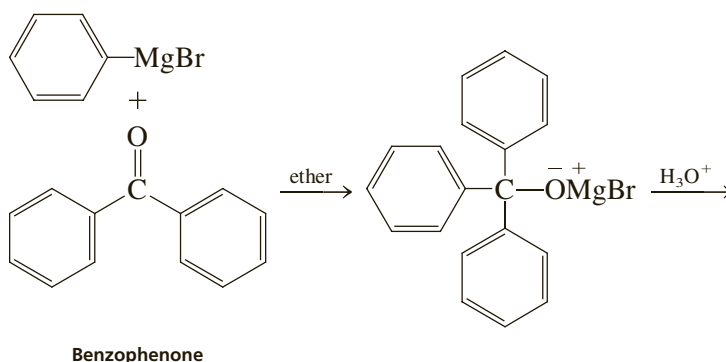


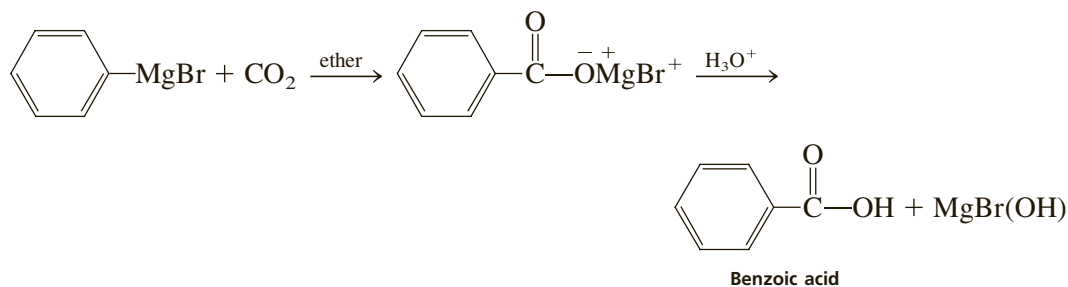**Benzoic Acid
(Experiment 35B)**

The alkyl portion of the Grignard reagent behaves as if it had the characteristics of a **carbanion**. We may write the structure of the reagent as a partially ionic compound:

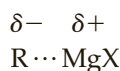


This partially bonded carbanion is a Lewis base. It reacts with strong acids, as you would expect, to give an alkane:

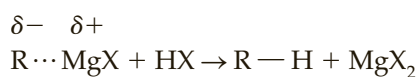


Any compound with a suitably acidic hydrogen will donate a proton to destroy the reagent. Water, alcohols, terminal acetylenes, phenols, and carboxylic acids are all acidic enough to bring about this reaction.


The Grignard reagent also functions as a good nucleophile in nucleophilic addition reactions of the carbonyl group. The carbonyl group has electrophilic character at its carbon atom (due to resonance), and a good nucleophile seeks out this center for addition.

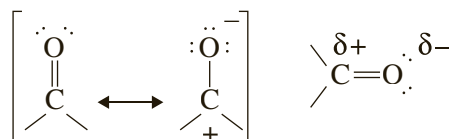

The magnesium salts produced form a complex with the addition product, an alkoxide salt. In a second step of the reaction, these must be hydrolyzed (protonated) by addition of dilute aqueous acid:

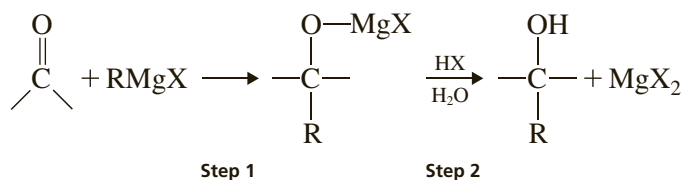


The Grignard reaction is used synthetically to prepare secondary alcohols from aldehydes and tertiary alcohols from ketones. The Grignard reagent will react with esters twice to give tertiary alcohols. Synthetically, it also can be allowed to react with carbon dioxide to give carboxylic acids and with oxygen to give hydroperoxides:

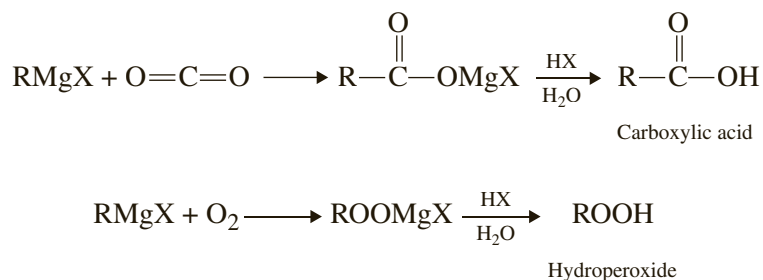

Because the Grignard reagent reacts with water, carbon dioxide, and oxygen, it must be protected from air and moisture when it is used. The apparatus in which the reaction is to be conducted must be kept scrupulously dry (recall that 18 mL of H₂O is 1 mole), and the solvent must be free of water, or anhydrous. During the reaction, the flask must be protected by a calcium chloride drying tube. Oxygen should also be excluded. In practice, this can be done by allowing the solvent ether to reflux. This blanket of solvent vapor keeps air from the surface of the reaction mixture.

In the experiment described here, the principal impurity is **biphenyl**, which is formed by a heat- or light-catalyzed coupling reaction of the Grignard reagent and unreacted bromobenzene. A high reaction temperature favors the formation of this product. Biphenyl is highly soluble in petroleum ether, and it is easily separated from triphenylmethanol. Biphenyl can be separated from benzoic acid by extraction.

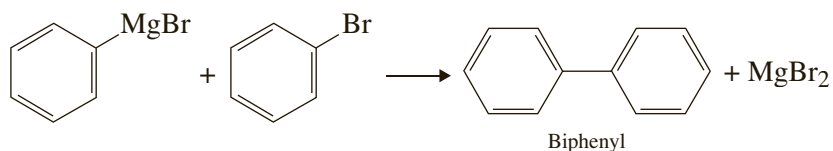

REQUIRED READING

<i>Review:</i> Technique 8	Section 8.3
Technique 11	Section 11.3
Technique 12	Sections 12.5, 12.7, 12.9, 12.11
Technique 25	Section 25.5

SPECIAL INSTRUCTIONS

This experiment must be conducted in one laboratory period either to the point after which benzophenone is added (Experiment 35A) or to the point after which the Grignard reagent is poured over dry ice (Experiment 35B). The Grignard reagent cannot be stored. This reaction involves the use of diethyl ether, which is extremely flammable. Be certain that no open flames of any sort are in your vicinity when you are using ether.

During this experiment, you will need to use *anhydrous* diethyl ether, which is usually contained in metal cans with a screw cap. You are instructed in the experiment to transfer a small portion of this solvent to an Erlenmeyer flask. Be certain to minimize exposure to atmospheric water. Always recap the container after use. Do not use solvent-grade ether because it may contain some water.

All students will prepare the Grignard reagent, phenylmagnesium bromide. At the option of the instructor, you should proceed to either Experiment 35A (triphenylmethanol) or Experiment 35B (benzoic acid).


SUGGESTED WASTE DISPOSAL

All aqueous solutions should be placed in the aqueous waste container. Be sure to decant these solutions away from any magnesium chips. The unreacted magnesium chips should be placed in a solid waste container.

Place all ether solutions in the container for nonhalogenated liquid wastes. Likewise, the mother liquor from the crystallization using isopropyl alcohol (Experiment 35A) should also be placed in the container for nonhalogenated liquid wastes.

PROCEDURE

Preparation of the Grignard Reagent: Phenylmagnesium Bromide

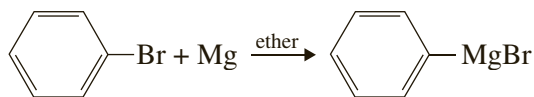

Preparation of Glassware


All glassware used in a Grignard reaction must be *scrupulously* dried. Surprisingly large amounts of water adhere to the walls of glassware, even glassware that is apparently dry. For this experiment, make sure all pieces of glassware have been rinsed with acetone and allowed to dry for at least two days. *If the equipment has been dried in this manner, then it is not necessary to dry the equipment in an oven.* Dry the following pieces of equipment before doing this experiment: a 20-mL round-bottom flask, two 5-mL conical vials, a 25-mL Erlenmeyer flask, a Claisen head, a syringe, and a calibrated Pasteur pipette (0.5-mL and 1.0-mL calibration marks) for use in dispensing ether. If, after drying as described, signs of water are still visible in the apparatus, dry the equipment in an oven. Prepare a drying tube with anhydrous calcium chloride.

Obtain 0.15 g of *shiny* magnesium turnings and place them in the *dry* round-bottom flask. Place a small *dry* magnetic stirring bar into the flask. Assemble the remainder of the apparatus, as shown in the figure. Seal off the open end of the Claisen head with a rubber septum.

Formation of the Grignard Reagent

Transfer about 10 mL of *anhydrous diethyl ether* into a dry 25-mL Erlenmeyer flask and stopper the flask. Use the flask to store your dry ether during the course of this experiment. During the experiment, remove the ether from this flask with a dry calibrated Pasteur pipette.

Place 0.70 mL of bromobenzene ($MW = 157.0$) into a preweighed 5-mL conical vial and determine the weight of the material transferred. Add 4.0 mL of *anhydrous* ether to the vial and mix the liquids. After the bromobenzene dissolves, withdraw about 0.8 mL of this solution into the syringe and cap the vial. You will need to save the remainder of the bromobenzene/ether solution for later use; recap the vial between uses. After inserting the syringe needles through the rubber septum, add 0.8 mL of the bromobenzene solution to the magnesium in the round-bottom flask. Position the apparatus just above the hot plate (about 60°C) and stir the mixture gently to avoid throwing the magnesium onto the side of the flask. You should begin to notice the evolution of bubbles, from the metal surface, which signals that the reaction is starting. It will probably be necessary to heat the mixture to start the reaction. Because ether has a low boiling point (35°C), it may be sufficient to heat the flask by placing it just above the hot plate. Check to see if the bubbling action continues after the apparatus is removed from the heat. The reaction should start, but if you experience difficulty, proceed to the next paragraph.

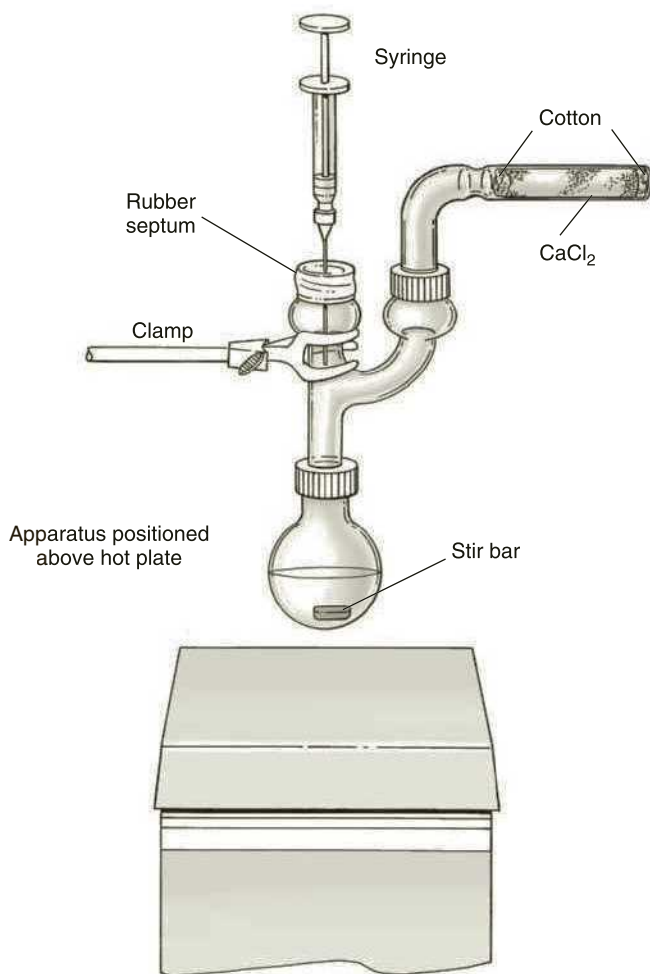

Apparatus for Experiment 35.

benzene solution and remove the flask from the hot plate. Ideally, the mixture will boil without the application of external heat. *If the reflux slows or stops, it is important that you heat the mixture.* As the reaction proceeds, you should observe the gradual disintegration of the magnesium metal. When all the bromobenzene has been added, place 2.0 mL of *anhydrous* ether in the vial that originally contained the bromobenzene solution, draw it into the syringe, and add the ether to the reaction mixture. Add more anhydrous ether to replace any that is lost during the reflux period. After a period of about 30 minutes from the beginning of the addition of bromobenzene, most or all of the magnesium should have reacted. Cool the mixture to room temperature. As your instructor designates, go to either Experiment 35A or Experiment 35B.

Optional Steps

You may need to employ one or more of the following procedures if heating fails to start the reaction. If you are experiencing difficulty, remove the syringe and rubber septum. Place a *dry* glass stirring rod into the flask and gently twist the stirring rod to crush the magnesium against the glass surface. Reattach the rubber septum and again heat the mixture. Repeat the crushing procedure several times, if necessary, to start the reaction. If the crushing procedure fails to start the reaction, then add one small crystal of iodine to the flask. Again, heat the mixture gently. The most drastic action, other than starting over again, is to prepare a small sample of the Grignard reagent in a test tube. When this reaction is started, it is added to the main reaction mixture in the flask.


Completing the Grignard Preparation

When the reaction has started, you should observe the formation of a brownish-gray, cloudy solution. Remove more of the bromobenzene/ether solution from the storage vial with the syringe and add the solution slowly over a period of 15 minutes. Refill the syringe as necessary until all the solution has been added to the magnesium metal. It may be necessary to heat the mixture occasionally with the hot plate during the addition, but if the reaction becomes too vigorous, slow the addition of the bromo-

35A EXPERIMENT 35A

Triphenylmethanol

PROCEDURE

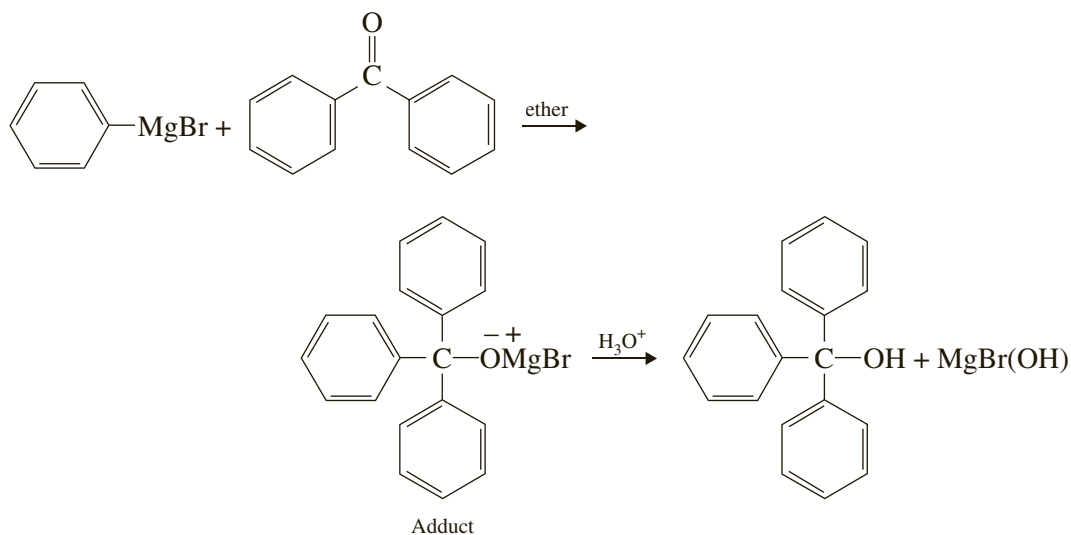**Addition of Benzophenone**

While the phenylmagnesium bromide solution is being heated and stirred under reflux, make a solution of 1.09 g benzophenone in 2 mL of *anhydrous* ether in a 5-mL conical vial. Cap the vial until the reflux period is over. Once the Grignard reagent is cooled to room temperature, draw some of the benzophenone solution into the syringe. Add this solution as rapidly (but not all at once) as possible to the stirred Grignard reagent. Do not add the solution so rapidly that the solution begins to boil. Add the remainder of the benzophenone solution with the syringe. Once the addition has been completed, cool the mixture to room temperature. The solution turns red and then gradually solidifies as the adduct is formed. When stirring is no longer effective, remove the syringe and septum and stir the mixture with a spatula. Rinse the vial that contained the benzophenone solution with about 1 mL of anhydrous ether and add it to the mixture. Remove the reaction flask from the apparatus and cap it. Occasionally stir the contents of this flask. Recap the flask when it is standing to avoid contact with water vapor. The adduct should be fully formed after about 15 minutes. You may stop here.

Hydrolysis

Add 6.0 mL of 6M hydrochloric acid (*dropwise at first*) to neutralize the reaction mixture. The acid converts the adduct to triphenylmethanol and inorganic compounds (MgX_2). Any unreacted magnesium will react with the acid to evolve hydrogen gas. Use a spatula to break up the solid while adding the hydrochloric acid. You may need to cap the flask and shake it vigorously to dissolve the solid. Because the neutralization procedure evolves heat, some ether will be lost due to evaporation. You should add enough additional ether to maintain at least a 10-mL volume in the upper organic phase. Eventually, you should obtain two distinct

layers: The upper ether layer will contain triphenylmethanol; the lower aqueous hydrochloric acid layer will contain the inorganic compounds. Make sure you have two distinct liquid layers, with no sign of any solid, before separating the layers. More ether or hydrochloric acid may be added, if necessary, to dissolve any remaining solid.

Transfer the entire contents of the reaction flask to a small separatory funnel, leaving the stirring bar behind. Use a small amount of ether to rinse the reaction flask and add this ether to the separatory funnel. If some solid material appears or if there are three layers, add more ether and hydrochloric acid to the separatory funnel and shake it. Continue adding small portions of ether and hydrochloric acid to the separatory funnel and shake it until everything dissolves. In some cases, it may be necessary to add more water instead of hydrochloric acid. Ultimately, you should have two distinct liquid layers with no sign of any solid, except possibly some magnesium. If a small amount of unreacted magnesium metal is present, you will observe bubbles of hydrogen being formed. You may remove the aqueous layer from the separatory funnel even though the magnesium is still producing hydrogen.

Separation and Drying

Drain the lower aqueous layer into a beaker. Pour the remaining ether layer that contains the triphenylmethanol product into a dry Erlenmeyer flask. Pour the aqueous layer back into the separatory funnel and reextract it with 5 mL of ether. Remove the lower aqueous phase and discard it. Combine the remaining ether phase with the first ether extract. Dry the ether solution with granular anhydrous sodium sulfate (Technique 12, Section 12.9).

Evaporation

Remove the dried ether solution from the drying agent by decanting it into a small Erlenmeyer flask and rinse the drying agent with more diethyl ether. Evaporate the solvent in a hood by heating the flask in a hot water bath at about 50°C (use an air or nitrogen stream to aid the evaporation process). After removal of the solvent, an oily solid should be left. This crude mixture contains the desired triphenylmethanol and the by-product, biphenyl. Most of the biphenyl can be removed by adding about 3 mL of *petroleum ether* (30 to 60°C). Petroleum ether is a mixture of hydrocarbons that easily dissolves the hydrocarbon biphenyl and leaves behind the alcohol triphenylmethanol. Do not confuse this solvent with diethyl ether (“ether”).


Heat the mixture slightly, stir it, and then cool the mixture to room temperature. Collect the triphenylmethanol by vacuum filtration on a Hirsch funnel and rinse it with small portions of petroleum ether (Technique 8, Section 8.3 and Figure 8.5). Air-dry the solid, weigh it, and calculate the percentage yield of the crude triphenylmethanol ($MW = 260.3$).

Crystallization

Crystallize your entire product from *hot* isopropyl alcohol in an Erlenmeyer flask using a hot plate as the heating source. Be sure to add the hot alcohol in small portions to the crude product. Add the hot solvent until the solid just dissolves. Then allow the flask to cool slowly. When it has cooled, place the flask in an ice bath to complete the crystallization. Collect the solid on a Hirsch funnel and wash it with a small amount of *cold* isopropyl alcohol. Set the crystals aside to air-dry. Report the melting point of the purified triphenylmethanol (literature value, 162°C) and recovered yield in grams. Submit the sample to the instructor in a properly labeled vial.

Spectroscopy

At the option of the instructor, determine the infrared spectrum of the purified material in a KBr pellet (Technique 25, Section 25.5). Your instructor may assign certain tests on the product you prepared. These tests are described in the Instructor's Manual.

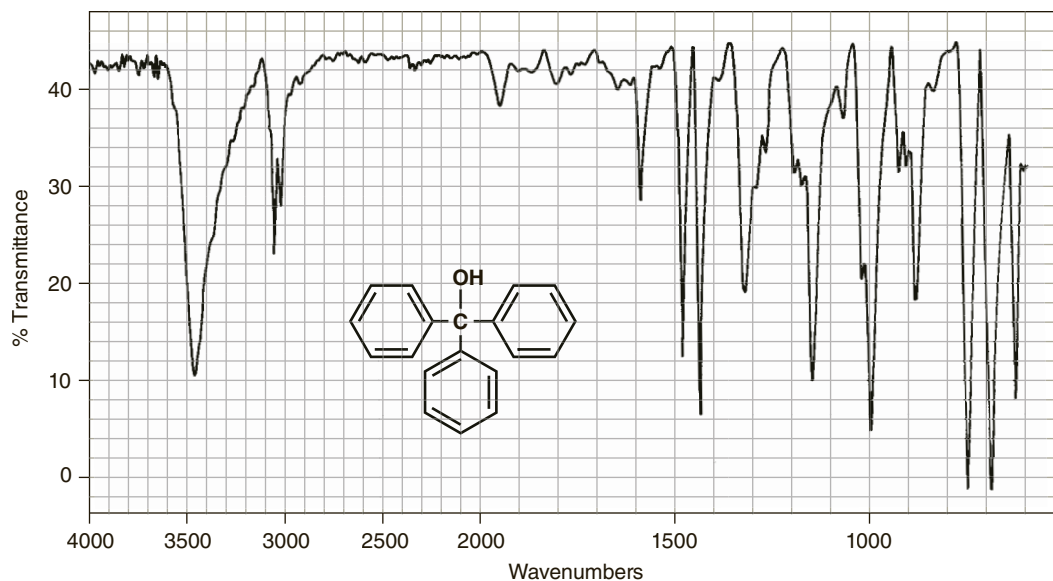


Infrared spectrum of triphenylmethanol, KBr.

35B EXPERIMENT 35B

Benzoic Acid

PROCEDURE

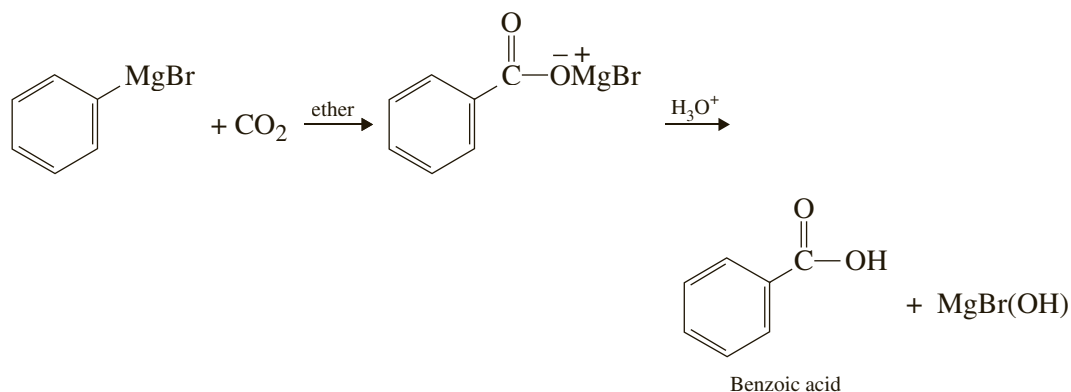

Addition of Dry Ice

When the phenylmagnesium bromide has cooled to room temperature, use a Pasteur pipette to transfer this reagent as quickly as possible to 4 g of crushed dry ice contained in a beaker. The dry ice should be weighed as quickly as possible to

avoid contact with atmospheric moisture. It need not be weighed precisely. Rinse the flask with 2 or 3 mL of anhydrous ether and add it to the beaker.

CAUTION



Exercise caution in handling dry ice. Contact with the skin can cause severe frostbite. Always use gloves or tongs. The dry ice is best crushed by wrapping large pieces in a clean, dry towel and striking them with a hammer or a wooden block. It should be used as soon as possible after crushing it to avoid contact with atmospheric water.

Cover the reaction mixture with a watch glass and let it stand until the excess dry ice has completely sublimed. The Grignard addition compound will appear as a viscous glassy mass.

Hydrolysis

Hydrolyze the Grignard adduct by slowly adding 10 mL of 6M hydrochloric acid, with stirring, to the beaker. Any remaining magnesium chips will react with acid to evolve hydrogen. If you have solid present (other than magnesium), try adding a little more ether. If the solid is insoluble in ether, try adding a little 6M hydrochloric acid solution. If neither seems to dissolve the solid, try adding a little water. Benzoic acid is soluble in ether, whereas the inorganic compounds (MgX_2) are soluble in the acid solution. Ultimately, you should have two distinct liquid layers in the beaker, with no sign of any solid, except possibly some magnesium. Transfer the liquid layers to a separatory funnel with a Pasteur pipette, leaving behind any residual magnesium.¹ If a separatory funnel is not available, you may use a centrifuge tube to separate the mixture. Add more ether to the beaker to rinse the beaker. Again, transfer the ether solution to the separatory funnel.

Isolation of the Product

Drain the lower aqueous layer and keep the upper ether layer in the separatory funnel. The aqueous phase contains inorganic salts and may be discarded. The ether layer contains the product, benzoic acid, and the by-product, biphenyl. Add 4 mL of 5% sodium hydroxide solution to the separatory funnel and shake it. Allow the layers to separate, drain the lower aqueous layer, and save this layer in a beaker. This extraction removes benzoic acid from the ether layer by converting it to the water-soluble sodium benzoate. The by-product, biphenyl, stays in the ether layer along with some remaining benzoic acid. Again, shake the remaining ether phase in the separatory funnel with a second 4-mL portion of 5% sodium hydroxide and drain the lower aqueous layer into the beaker with the first extract. Repeat the extraction process with a third portion (4 mL) of 5% sodium hydroxide and save the aqueous layer, as before. Discard the ether layer that contains the biphenyl impurity.

Heat the combined basic extracts with stirring on a hot plate (hot enough to boil the aqueous mixture) for about 5 minutes to remove any ether that may be dissolved in this aqueous phase. Stir the mixture as it is being heated. Ether is soluble in water to the extent of 7%. During this heating period, you may observe slight

¹ If it is necessary to set this experiment aside overnight, do not transfer the solution to the separatory funnel. Instead, transfer the solution to a 25-mL Erlenmeyer flask. Stopper the flask tightly. When you resume the experiment, transfer this solution to the separatory funnel using about 4 mL of ether to aid the transfer and proceed as instructed.

bubbling, but the volume of liquid will not decrease substantially. Unless the ether is removed before the benzoic acid is precipitated, the product may appear as a waxy solid instead of crystals.


Cool the alkaline solution and precipitate the benzoic acid by adding 5 mL of 6M hydrochloric acid with stirring. Cool the mixture in an ice bath. Collect the solid by vacuum filtration on a Hirsch funnel (Technique 8, Section 8.3, and Figure 8.5). The transfer may be aided and the solid washed with several small portions of cold water (total volume, 4 mL). Allow the crystals to dry thoroughly at room temperature at least overnight in an open container. Weigh the solid and calculate the percentage yield of benzoic acid ($MW = 122.1$).

Crystallization

Crystallize your entire product from *hot* water in an Erlenmeyer flask using a hot plate as the heating source. Be sure to add the hot water in small portions to the crude product. Add the hot water until the solid just dissolves. Then allow the flask to cool slowly. After the flask cools to room temperature, place it in an ice bath to complete the crystallization. Collect the solid on a small Hirsch funnel and wash it with a small amount of *cold* water. Set the crystals aside to air-dry at room temperature until the next laboratory period before determining the melting point of the purified benzoic acid (literature value, 122°C). Also determine the recovered yield in grams. Submit the sample to the instructor in a properly labeled vial.

Spectroscopy




At the option of the instructor, determine the infrared spectrum of the purified material in a KBr pellet (Technique 25, Section 25.5). Your instructor may assign certain tests on the product you prepared. These tests are described in the Instructor's Manual.

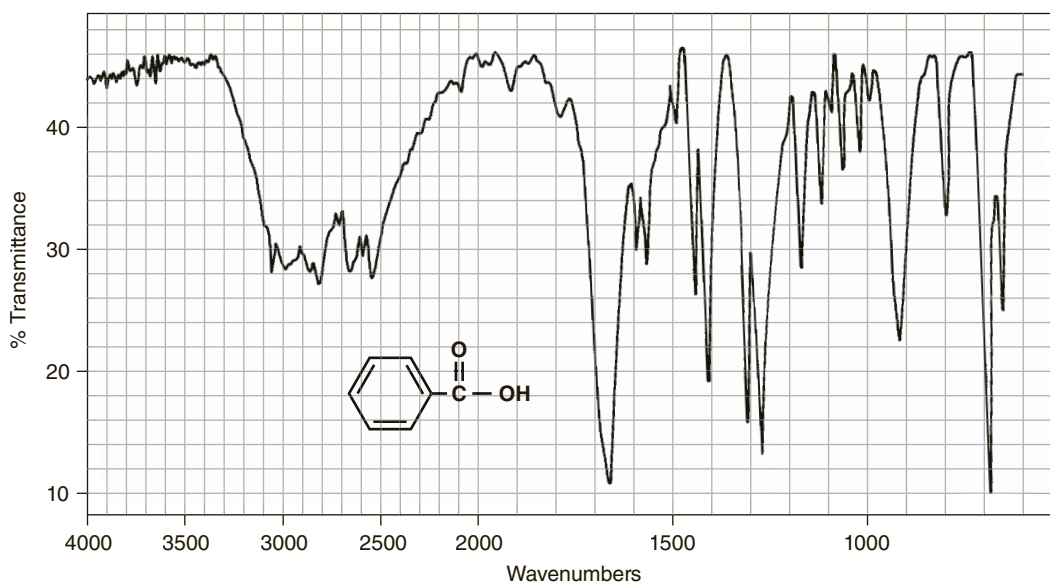

Infrared spectrum of benzoic acid, KBr.

QUESTIONS

1. Benzene is often produced as a side product during Grignard reactions using phenylmagnesium bromide. How can its formation be explained? Give a balanced equation for its formation.
2. Write a balanced equation for the reaction of benzoic acid with hydroxide ion. Why is it necessary to extract the ether layer with sodium hydroxide?
3. Interpret the principal peaks in the infrared spectrum of either triphenylmethanol or benzoic acid, depending on the procedure used in this experiment.
4. Outline a separation scheme for isolating either triphenylmethanol or benzoic acid from the reaction mixture, depending on the procedure used in this experiment.
5. Provide methods for preparing the following compounds by the Grignard method:

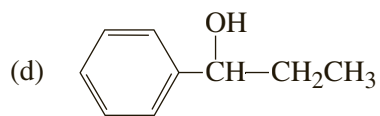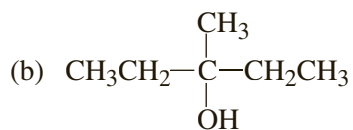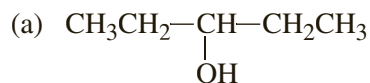

Aqueous-Based Organozinc Reactions

Organometallic reactions

Green chemistry


Extractions

Use of a separatory funnel

Gas chromatography

Spectroscopy


One of the most important categories of reactions in organic synthesis is the class of reactions that result in the formation of a carbon-carbon bond. Among these, one of the best-known reactions is the Grignard reaction, where an organomagnesium reagent is formed from an alkyl halide and then allowed to react with a variety of substances to form new molecules. The nucleophilic nature of the organomagnesium reagent is used in the formation of new carbon-carbon bonds. The equation shown illustrates this type of synthesis. The Grignard reaction is introduced in Experiment 35.

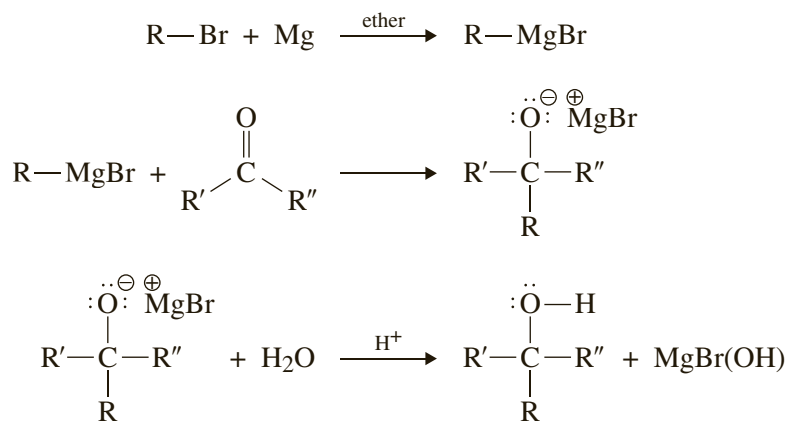

Because the organomagnesium reagent reacts with water, carbon dioxide, and oxygen, it must be protected from air and moisture when it is used. The apparatus in which the reaction is to be conducted must be scrupulously dry, and the solvent must be completely anhydrous. In addition, diethyl ether is required as a solvent; without the presence of an ether, the organomagnesium reagent will not form.

This experiment presents a variation on the basic idea of a Grignard synthesis, but one that does not use magnesium and that can be conducted in a mixed organic-aqueous solution. The reaction presented in this experiment is a variation on the Barbier-Grignard reaction, where zinc is used as the metal. A small amount of an ether, in this case tetrahydrofuran (THF), is still required for this reaction, but the principal component of the solvent system is water. By eliminating much of the organic solvent, this method can be used to illustrate some of the principles of

“green chemistry,” in which reactions are conducted under conditions that are less harmful to the environment than traditional chemical methods.

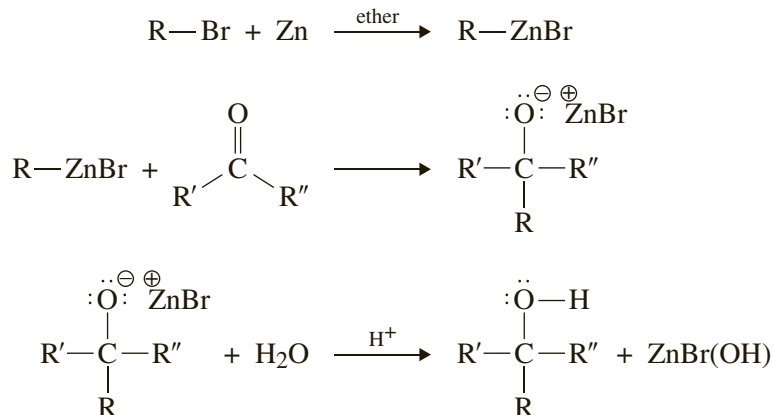

Although this organozinc method of synthesis is very similar to the Grignard reaction, there are also some interesting differences. The organozinc reagent is much more selective than the organomagnesium reagent, and rearrangements of the alkyl group attached to the metal are also possible. Whereas the formation of Grignard reagents from allylic halides is notoriously difficult, the formation of organozinc reagents seems to require that one begin with an allylic halide. A comparison of the structure of the products of this reaction with the structure of the starting alkyl halide can reveal some of this interesting chemistry.

REQUIRED READING

Review: Technique 8 Section 8.3
 Technique 7 Section 7.10
 Technique 12 Sections 12.5, 12.8, 12.9, 12.11
 Technique 22
 Technique 25 Sections 25.2, 25.4
 Technique 26 Section 26.1
 Technique 27 Section 27.1

SPECIAL INSTRUCTIONS

This reaction involves the use of allyl bromide, a substance that is volatile and may also be a **lachrymator**. Be certain to dispense this material under the hood. Do not attempt to weigh this substance; determine the approximate volume of allyl bromide needed using the specific gravity provided in this experiment, and dispense the allyl bromide by volume using a calibrated pipette. Students should work in pairs for this experiment.

SUGGESTED WASTE DISPOSAL

All aqueous solutions should be placed in a waste container designated for the disposal of aqueous wastes.

PROCEDURE

ACTIVATED ZINC

Carefully weigh 1.31 g (0.02 moles) of zinc powder and add it to a small (10-mL) Erlenmeyer flask or beaker. Add 1 mL of 5% aqueous hydrochloric acid and allow the mixture to stand for 1 to 2 minutes. There will be a noticeable evolution of hydrogen gas during this time. At the end of this period, pour the entire mixture into a Hirsch funnel and isolate the zinc by vacuum filtration. Rinse the zinc with 1 mL of water, followed by 1 mL of ethanol and 1 mL of diethyl ether. The zinc should be ready to use for the procedure, as described below.

PREPARATION AND REACTION OF THE ORGANOZINC REAGENT

Add 10 mL of saturated aqueous ammonium chloride solution to a 25-mL round-bottom flask. Add 1.31 g zinc powder (0.02 moles) and a stirring bar to the flask. Attach an air condenser to the flask and begin continuous stirring while adding the remaining reagents. Carefully weigh 0.86 g (0.01 moles) of the 3-pentanone. Add the ketone and 1.6 mL of tetrahydrofuran to a test tube and add this solution dropwise to the zinc/ NH_4Cl solution. The rate of addition should be about one drop per second. Note that this addition can be made by dropping the solution carefully down the opening in the air condenser; use a Pasteur pipette to add the solution. Allow the solution to stir for 10 to 15 minutes, giving time for the carbonyl compound to form a complex with the zinc. Add 2.4 g (0.02 moles—use the specific gravity 1.398 g/mL to estimate the volume required) of allyl bromide (3-bromopropene) to the stirring solution. *Be sure to dispense this reagent in the hood!* The rate of addition should be about one drop per second. Add the halide carefully by dropping it down the opening in the air condenser. Allow the reaction mixture to stir for 1 hour.

Set up a vacuum filtration apparatus with a Hirsch funnel. Decant the liquid from the reaction mixture through the Hirsch funnel. Rinse the round-bottom flask with approximately 1 mL of diethyl ether and pour the liquid into the Hirsch funnel. Using a second 1-mL portion of diethyl ether, rinse the solid that has collected in the Hirsch funnel. Discard the solid. Prepare a filter-tip pipette and transfer the liquid that was collected in the vacuum filtration into a separatory funnel. Use 1 mL of diethyl ether to rinse the inside of the filter flask and use the filter-tip pipette to transfer this liquid to the separatory funnel. Shake the separatory funnel gently to extract the organic material from the aqueous layer to the ether layer. Drain the lower (aqueous) layer into a 50-mL Erlenmeyer flask. *Do not discard this aqueous layer.* Collect the upper (organic) layer from the separatory funnel into a 25-mL Erlenmeyer flask (remember to collect the upper layer by pouring it from the top of the separatory funnel). Replace the aqueous layer in the separatory funnel and wash it with a 2-mL portion of ether. Separate the layers, save the aqueous layer in the same 50-mL Erlenmeyer flask as before, and combine the ether layer with the ether solution collected in the previous extraction. Repeat this extraction process of the aqueous phase one more time using a fresh 2-mL portion of ether. Dry the combined ether extracts with 3–4 microspatulafuls of anhydrous sodium sulfate (see Technique 12, Section 12.9). Stopper the Erlenmeyer flask with a cork and allow it to stand for at least 15 minutes (or overnight).

Use a filter-tip pipette to transfer the dried liquid to a clean, preweighed Erlenmeyer flask. Use a small amount of ether to rinse the inside of the original flask and

add this ether to the dried liquid. Evaporate the ether with a rotary evaporator or under a gentle stream of air. When the ether has evaporated completely, reweigh the flask to determine the yield of product. If it should be necessary to store your final product, use Parafilm[®] to seal the container.

Prepare a sample of your final product for analysis by gas chromatography. Determine the infrared spectrum and both proton and ¹³C NMR spectrum of your product. Use these spectra to determine the structure of your product. In your laboratory report, include an interpretation of each spectrum, identifying the principal absorption bands and demonstrating how the spectrum corresponds to the structure of your compound. Submit your sample in a labeled vial with your laboratory report.

QUESTIONS

1. Write *balanced* chemical equations for the formation of a substance that you prepared in this experiment.
2. Outline a series of chemical equations to show how your product could have been prepared using a Grignard reaction. Be sure to show the structures of all starting materials and intermediates.
3. Draw the structure of the product that would have been formed if benzaldehyde had been used in place of 3-pentanone in this experiment.
4. When benzaldehyde is used as the carbonyl compound in this experiment, the CH₂ peak in the proton NMR spectrum appears as *two* separate, complex resonances. Explain why this is observed. (Hint: see Technique 26, Section 26.16.)


Sonogashira Coupling of Iodosubstituted Aromatic Compounds with Alkynes using a Palladium Catalyst

Green chemistry


Organometallic chemistry

Palladium-catalyzed reaction


In this experiment, we will conduct some modern organic chemistry using a palladium catalyst. It is a rare opportunity for students in undergraduate laboratories to experience this powerful chemistry. We will react the iodosubstituted aromatic compounds, shown below, with 1-pentyne, 1-hexyne, or 1-heptyne in the presence of the catalysts, palladium acetate and cuprous iodide, to yield 4-substituted-1-pentynyl, 4-substituted-1-hexynyl, or 4-substituted-1-heptynylaromatic compounds. This reaction is called the Sonogashira coupling reaction.¹ The reaction will be carried out in refluxing 95% ethanol as the solvent. In addition, piperazine will be employed both as a base and as a hydride donor.

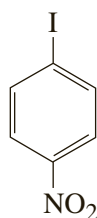


1-iodo-4-nitrobenzene

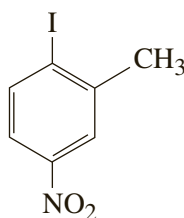

2-iodo-5-nitrotoluene

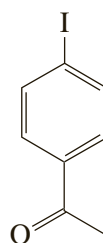

4-iodoacetophenone

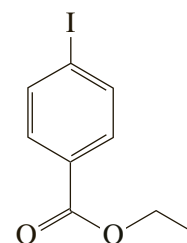

ethyl 4-iodobenzoate

Background

Palladium-catalyzed reactions can be used to connect the terminal end of an alkyne and aromatic iodide, as shown in the reaction below.² They are useful in industry and are widely employed in the academic arena. The experiment presented here was adapted from an article by Goodwin, Hurst, and Ross.³ The mechanism shown is for the coupling of 1-iodo-4-nitrobenzene with 1-pentyne. Small amounts of a dimer obtained from the coupling of the 1-alkynes are also formed in these


¹ a) Takahashi, S., Kuroyama, Y., Sonogashira, K., Hagihara, N. *Synthesis*, **1980**, 627–630.

b) Thorand, S., Krause, N. *J. Org. Chem.*, **1998**, 63, 8551–8553.

² Brisbois, R. G., Batterman, W. G., Kragerud, S. R. *J. Chem. Ed.* **1997**, 74, 832–833.


³ Goodwin, T. E., Hurst, E. M., Ross, A. S. *J. Chem. Ed.* **1999**, 76, 74–75. Experiment developed by Brogan, H., Engles, C., Hanson, H., Phillips, S., Rumberger, S., Lampman, G. M., Western Washington University, Bellingham, WA.

reactions. It is likely that the dimers result from the formation of the copper intermediate (Step 3 of the mechanism). Thus, reactions involving 1-pentyne yield some 4,6-decadiyne.

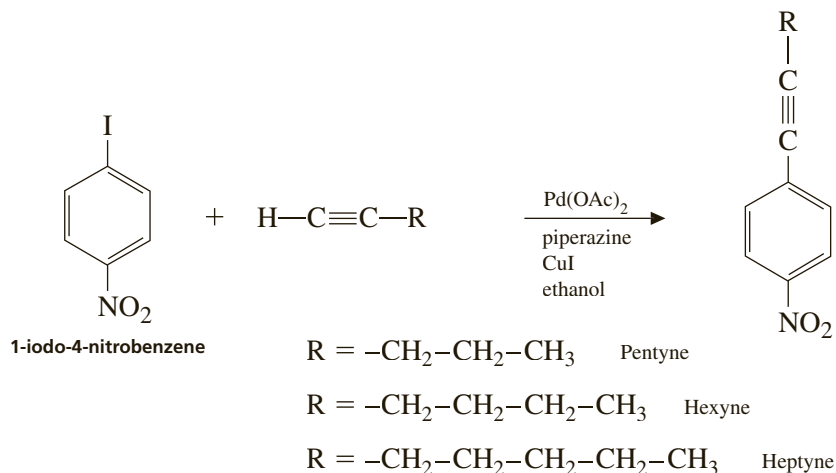


The mechanism is thought to proceed in six steps, as shown below.


Step 1: Transfer of hydride from piperazine to palladium

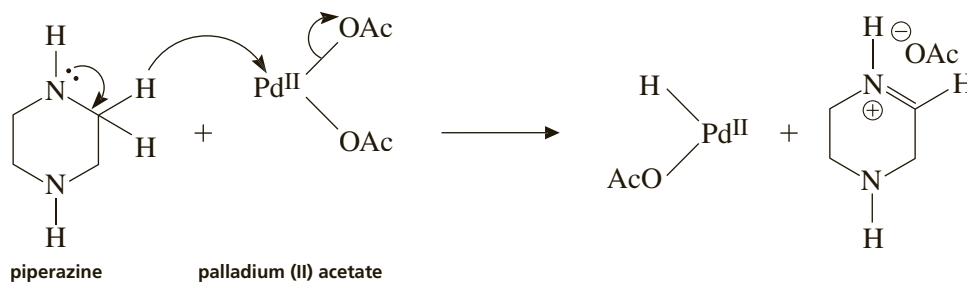


Step 2: Reduction of Pd(II) to Pd⁰ by removal of HOAc with piperazine

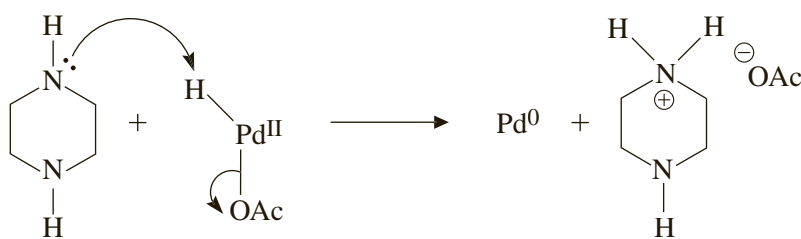


The Pd⁰ is probably complexed with piperazine ligands (L).

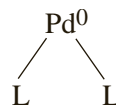



Step 3: Preparation of cuprate

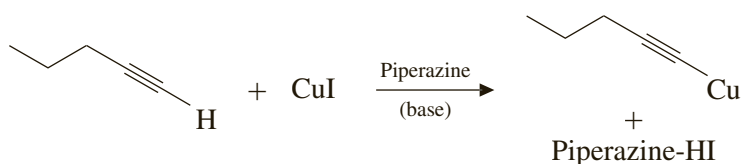

Step 4: Oxidative addition

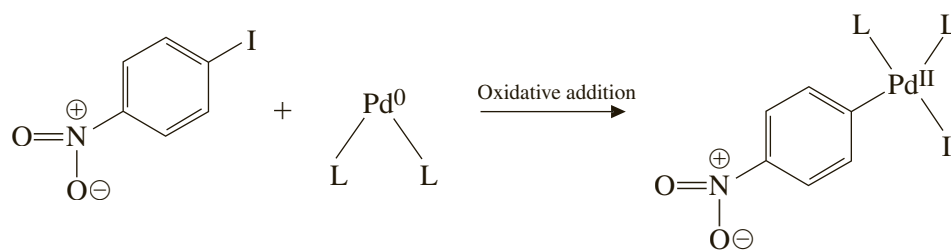

Step 5: Coupling of cuprate to the palladium complex

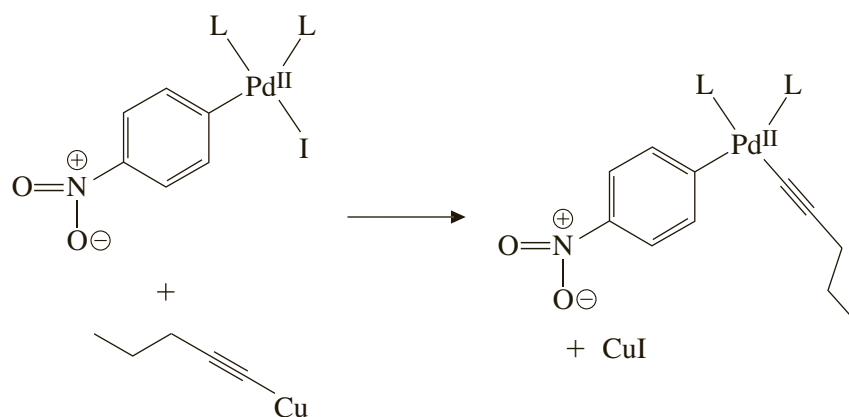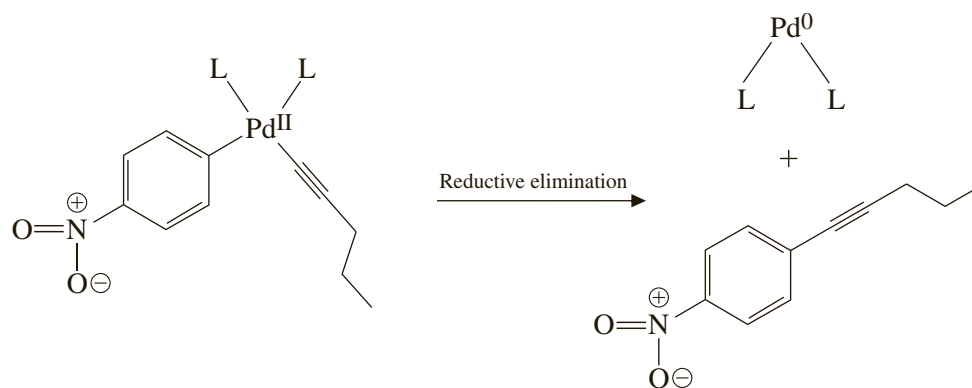Step 6: Reductive elimination forms the product and regenerates Pd⁰**REQUIRED READING**

Review: Techniques 5, 6, 7, 12, 19, 25, and 26

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous wastes in the container for aqueous waste. Place the organic waste in the nonhalogenated organic waste container. Place the halogenated waste in the appropriate waste container.

NOTES TO THE INSTRUCTOR

It is suggested that students work in pairs for this experiment. A rotary evaporator and vacuum pump are useful items to have in the laboratory, but if they are unavailable, the solvents can be removed easily by heating the sample, while blowing air or nitrogen on the sample.

The Sonogashira reaction works best when electron withdrawing functional groups are attached to the aromatic ring. Thus, the four compounds shown above work well employing a 30-minute reaction period. These compounds contain nitro, acetyl, and carboethoxy functional groups, along with the iodo group. When electron-releasing groups such as methoxy are attached to the ring, the reaction is much slower and requires much longer reaction times. We have found success with the less reactive compounds employing microwave technology. If your laboratory includes a commercial microwave reactor, such as the CEM Explorer, you can achieve excellent success with 4-iodoanisole (1-iodo-4-methoxybenzene) using the optional procedure.

PROCEDURE

Preparation of the Reaction Mixture

Add 0.200 mmol of one of the four iodo substrates shown to a 25-mL round-bottom flask. Use a 4-place analytical balance for weighing the substrates and all of the materials listed next. Now add 55 mg of piperazine and a clean magnetic stir bar to the flask. Add 1.25 ml of 95% ethanol to the flask to dissolve the materials. Now add 16.5 mg of palladium (II) acetate and 10 mg of copper (I) iodide to the flask. Finally, use an automatic pipette to dispense 70 μ L of 1-pentyne, 1-hexyne, or 1-heptyne, depending on which alkyne you were assigned, to the round-bottom flask. Attach a water-cooled condenser to the flask. Heat the contents at reflux for 30 minutes on a hot plate, with stirring.

After the solution has been refluxed for 30 minutes, allow the mixture to cool for a few minutes. Remove the flask and remove the ethanol on a rotary evaporator.⁴ When using the rotary evaporator, be sure to spin the flask rapidly and don't heat the water in the water bath. There may be a tendency for the sample to "bump." When it appears that the ethanol has been removed, attach the flask to a vacuum pump for at least 3 minutes to remove the remaining ethanol and any dimer formed in the reaction. When the ethanol has been successfully removed, add 1 mL of methylene chloride to the flask followed by 0.2 g of silica gel. Swirl the flask to ensure that most of the liquid is adsorbed onto the silica gel. Put the flask back onto the rotary evaporator and remove the methylene chloride⁴. Your product is now adsorbed onto the silica, yielding a dry, free-flowing solid. Use a spatula to break up the silica containing your product. Pour the solid onto a piece of paper and keep it handy until you have made up the column.

Column Chromatography

Prepare a silica gel column for chromatography using a 10-mL Pyrex disposable cleanup/drying column (Corning #214210 available from Fisher #05-722-13; the column is about 30 cm long and 1 cm in diameter). Push some cotton down into the bottom using a glass rod, but don't pack the cotton too tightly. The cotton must be

⁴ An alternative procedure for removing the ethanol and CH_2Cl_2 solvent, is to blow air on the sample. Allow at least 10 to 15 minutes at 50°C for removal of the ethanol.

tight enough to keep the silica gel from leaking out of the bottom of the column, but not too tight or it will reduce the flow of solvent. Add silica gel⁵ until it is about 5 cm from the top of the column.

Construct a funnel from a disposable plastic Pasteur pipette in order to add the sample to the top of the chromatography column. To make the funnel, first cut off the top of a 1-mL plastic pipette and also remove most of the tip to make a small funnel (your instructor will demonstrate this). Pour the silica sample containing your adsorbed product from the weighing paper into the top of the silica gel column using your funnel. The solid now resides at the top of your chromatography column. Obtain 10 mL of hexanes and 20 mL of CH₂Cl₂. First pass the 10 mL of hexane through the column, in portions, to wet the silica, and collect the eluent in a preweighed 100-mL round-bottom flask (obtain the flask from your instructor, and use a 4-place balance). Then pass the CH₂Cl₂ solvent through the column in portions while collecting the eluent into the same 100-mL flask. The column removes the palladium catalyst, which remains as a black substance at the top of the chromatography column.

Isolation of the Product

After all of the elutants have been collected in the round-bottom flask, attach the flask to the rotary evaporator and remove the solvent, under vacuum.⁴ (Be careful that the solvent doesn't bump up into the trap!) After removing the hexanes and CH₂Cl₂, attach the flask to a vacuum pump⁶ for about 3 minutes to ensure that all of the solvent and dimer⁷ have been removed from the product. Remove the flask and weigh it on the 4-place balance to determine the amount of product obtained. Calculate the percentage yield.

Analysis of the Product


Determine the NMR spectrum of the sample remaining in the 100-mL flask in CDCl₃. Add a few drops of CDCl₃ directly to the flask. Transfer the solution to the NMR tube with a Pasteur pipette. Put more drops of CDCl₃ into the flask, and transfer this to the NMR tube. Repeat until you are fairly certain that you have transferred most of your product to the NMR tube. Finally, if necessary, add enough CDCl₃ to bring the total height to about 50 mm. Run the NMR spectrum and interpret the patterns. Four reference spectra are shown in Figures 1, 2, 3, and 4. Figure 1 shows the spectrum for the product obtained from 1-iodo-4-nitrobenzene and 1-hexyne. Notice that the spectrum shows a triplet at 0.96 ppm, a sextet at 1.50 ppm, a quintet at 1.60 ppm, another triplet at 2.45 ppm, and 2 doublets—one at 7.50 ppm and one at 8.15 ppm. A trace of 5,7-dodecadiyne is observed at about 0.9, 1.4, and 2.2 ppm in the NMR spectrum. Be on the alert for a sharp singlet that may appear near 7.25 ppm for traces of chloroform (CHCl₃) present in the CDCl₃ solvent. Other NMR spectra are shown in Figures 2, 3, and 4. Compare your results to those shown in the figures when making the assignments for your sample.

The plan is to run the proton NMR spectrum, and then use your sample to obtain the infrared spectrum. Pour the contents of the NMR tube into a small test tube. Transfer a small amount of the CDCl₃ solution to a salt plate using a Pasteur pipette, blow on the plate to evaporate the solvent, and then determine the infrared


⁵ Fisher Chromatographic Silica Gel, 60-200 mesh, #S818-1, Davisil® Grade 62, type 150A°.

⁶ The vacuum pump is helpful to remove all traces of hexane and methylene chloride. In the NMR spectrum, these peaks appear at 0.9 ppm (triplet) and 1.3 ppm (multiplet). Any remaining CH₂Cl₂ appears at about 5.3 ppm (singlet). See notes to the instructor.


⁷ The vacuum pump helps remove any dimer present in the sample. Be sure to use a good quality vacuum pump to remove the dimer from the product. See notes to the instructor.

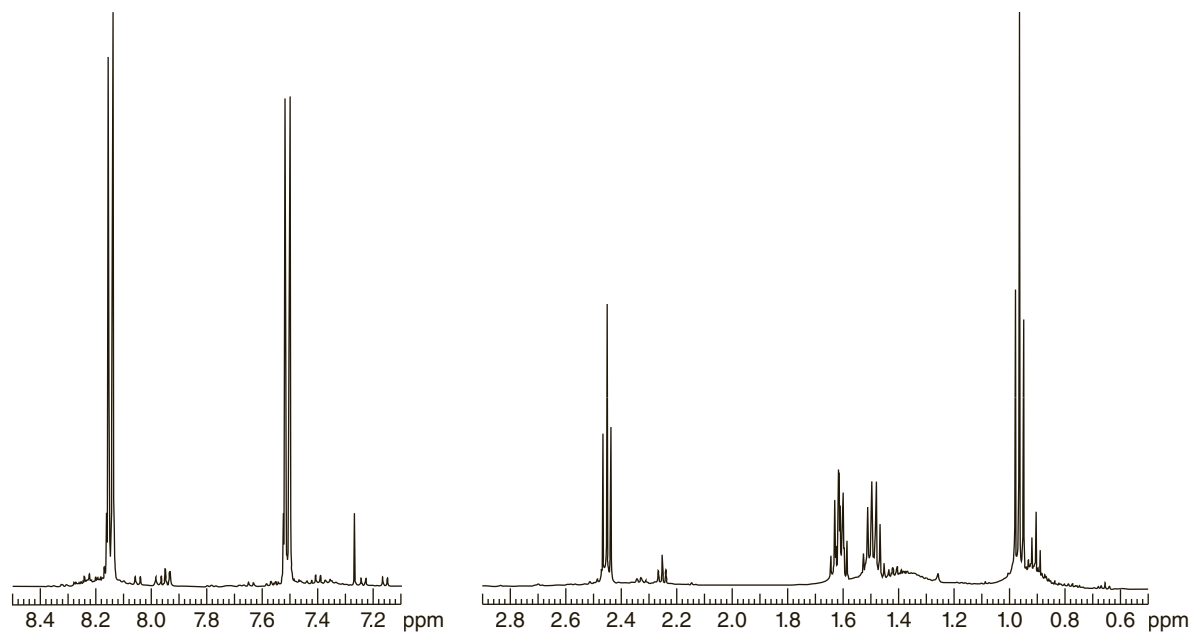**Figure 1**


500 MHz NMR spectrum of the product of 1-iodo-4-nitrobenzene and 1-hexyne. A trace of a dimer, 5,7-dodecadiyne, formed from 1-hexyne is observed at about 0.9 ppm (3H), 1.4 ppm (4H), and 2.2 ppm (2H) in the NMR spectrum. Traces of other impurities are also found in the spectrum. CHCl₃ appears at about 7.25 ppm.

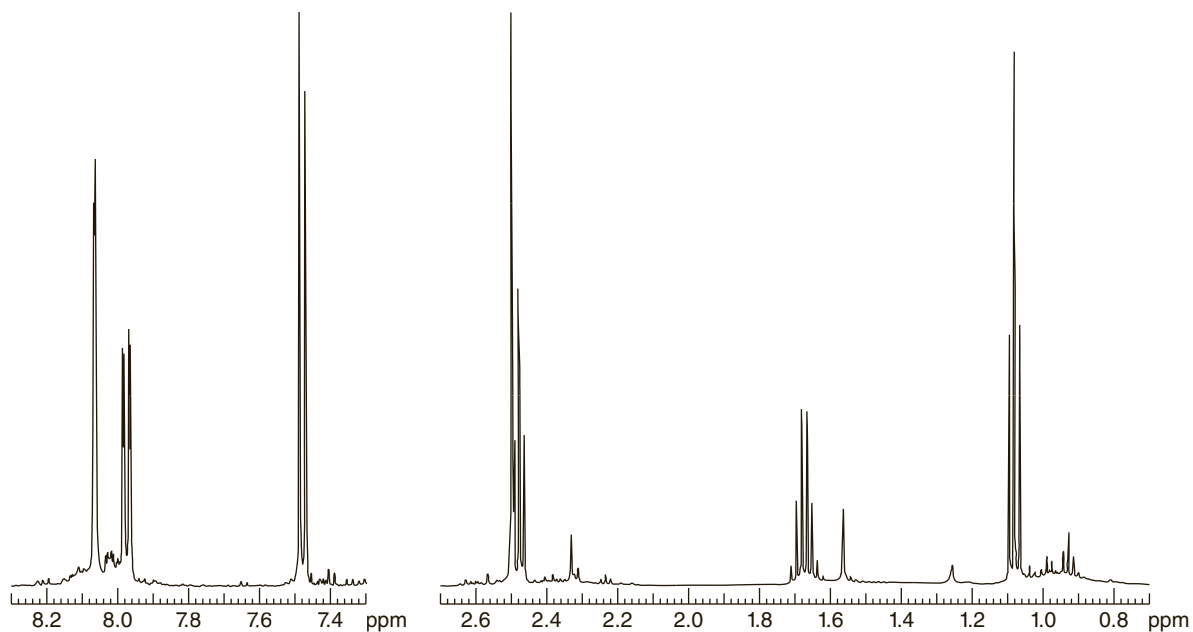**Figure 2**


500 MHz NMR spectrum of the product of 1-iodo-2-methyl-4-nitrobenzene and 1-pentyne. Notice that the singlet for the methyl group partially overlaps with the triplet at 2.5 ppm.

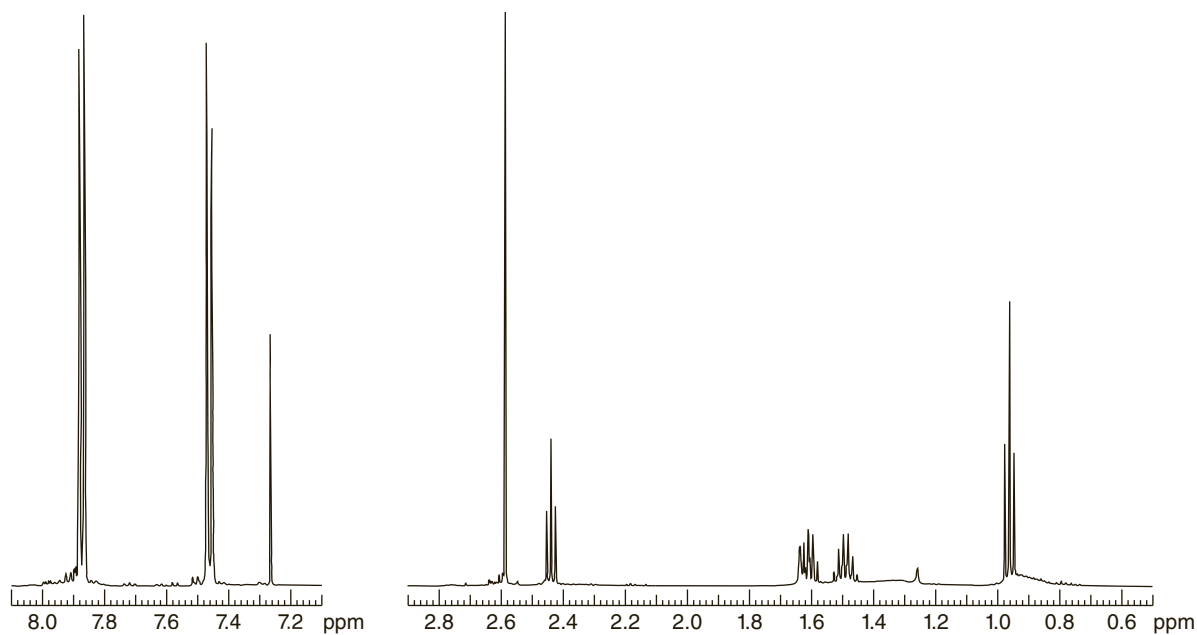**Figure 3**

500 MHz NMR spectrum of the product of 4-iodoacetophenone and 1-hexyne. CHCl_3 appears at about 7.25 ppm.

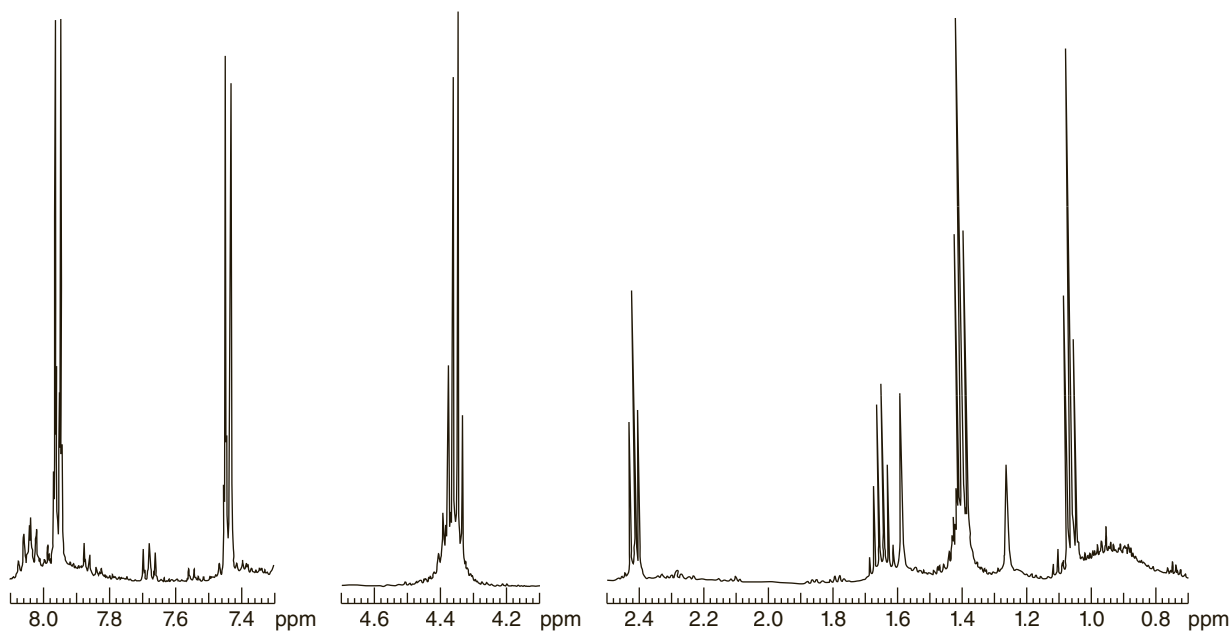**Figure 4**

500 MHz NMR spectrum of the product of ethyl 4-iodobenzoate and 1-pentyne. The $-\text{CH}_2-$ in the ethyl group appears as a quartet at 4.4 ppm, while the CH_3 group in the ethyl group appears as a triplet at 1.4 ppm. The triplet at 1.05 ppm, sextet at 1.65 ppm, and triplet at 2.40 ppm are assigned to the $-\text{CH}_2-\text{CH}_2-\text{CH}_3$ chain. The pair of doublets at 7.45 and 7.95 ppm are assigned to the *para*-disubstituted benzene ring. Impurity peaks appear at 0.95 (broad) ppm, 1.25 ppm, and 1.60 ppm and along with some miscellaneous small impurities appearing in the aromatic ring region.

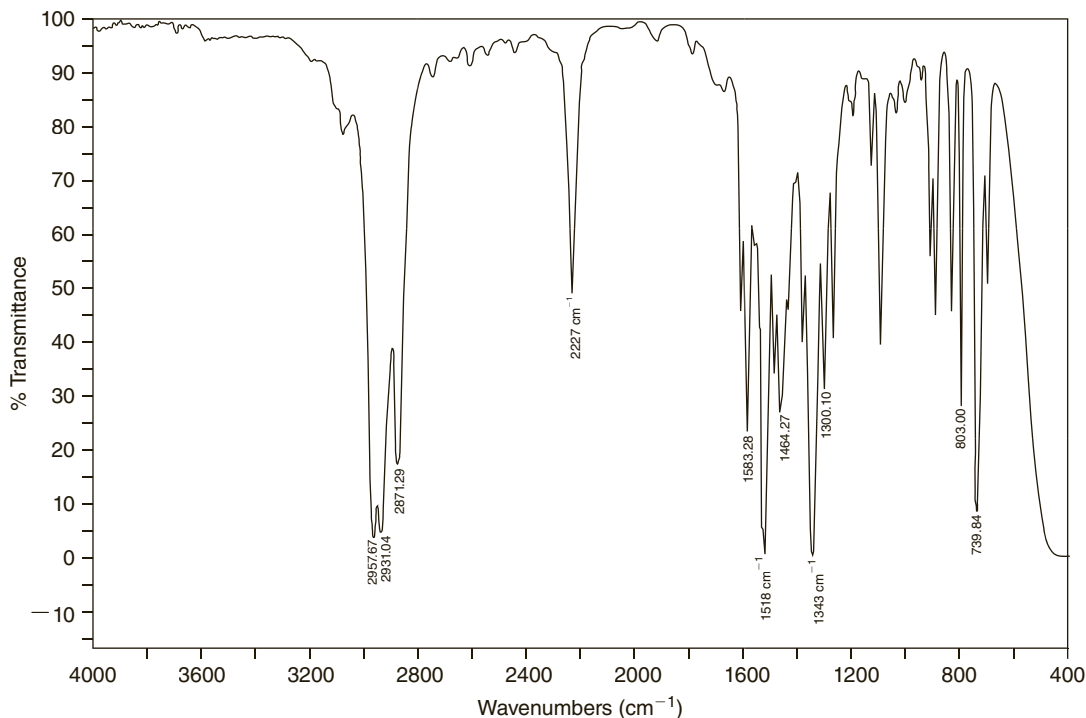**Figure 5**

Infrared spectrum of the product of 1-iodo-2-methyl-4-nitrobenzene and 1-hexyne. The sharp peak at 2227 cm^{-1} is assigned to the triple bond in 1-hexynyl-2-methyl-4-nitrobenzene, and the two sharp peaks at 1518 and 1343 cm^{-1} are assigned to the nitro group.

spectrum. Make sure that the CDCl_3 has evaporated before determining the infrared spectrum. The infrared spectrum for the product obtained from 1-iodo-2-methyl-4-nitrobenzene and 1-hexyne is shown in Figure 5 for comparison. A sharp peak at about 2227 cm^{-1} is observed for the triple bond, as well as two sharp peaks at 1518 and 1343 cm^{-1} for the NO_2 group. Assign peaks for your compound.

OPTIONAL PROCEDURE USING MICROWAVE TECHNOLOGY⁸

Reaction⁹

Add 0.0573 g (0.24 mmol) of 4-iodoanisole, 0.0120 g of palladium black powder, 0.1460 g of 40% potassium fluoride on alumina (Aldrich Chemical Co. #316385), 0.0317 g triphenylphosphine, 0.0410 g of cuprous iodide, 1 mL of 95% ethanol, and 70 μL of 1-pentyne to a standard microwave tube (12 mL). Add a stir bar recommended by the manufacturers of the microwave reactor. Cap the microwave tube securely with one of the caps supplied by the manufacturer of the microwave unit.

Microwave Instrument Conditions

Using the software supplied by the manufacturer, set the instrument to run at 100°C for 30 min with stirring on high.

⁸ Microwave apparatus: CEM Explorer, CEM Corp, 3100 Smith Farm Road, Mathews, NC 28106-0200.

⁹ Kabalka, G. W., Wang, L., Namboodiri, V., Pagni, R. M. Rapid microwave-enhanced, solventless Sonogashira coupling reaction on alumina, *Tetrahedron Letters*, **2000**, *41*, 5151–5154.

Workup Procedure

Following the 30 min reaction period and cooling period, add another 1-mL portion of 95% ethanol and vacuum-filter the mixture (see Technique 8, Figure 8.5) through a Hirsch funnel with filter paper to remove all of the solids present in the reaction tube. Aid the transfer process by using about 3 mL of 95% ethanol

Purification Procedure



Using a 1-mL pipette, transfer the liquid contents in the filter flask to a preweighed 25-mL round-bottom flask. Remove the ethanol, under vacuum, with a rotary evaporator. When it appears that the ethanol has been removed on the rotary evaporator, remove the flask and attach the flask to a good vacuum pump source to remove the remaining ethanol and any dimer (4,6-decadiyne) that may have formed in the reaction from the 1-pentyne. Continue pumping on the flask for at least 3 minutes. Release the vacuum, remove the flask, and reweigh the flask to determine the amount of product obtained. Calculate the theoretical yield for the reaction.


NMR Spectroscopy

Add about 0.7 mL of CDCl_3 to the sample in the flask. In most cases, you will find a small amount of undesired solid present that does not dissolve in the CDCl_3 . Prepare a filtering pipette (see Technique 8, Section 8.1C), draw up the CDCl_3 solution with a Pasteur pipette, and add it to the filtering pipette. Collect the solution in a small test tube. This filtering process should remove all or most of the solid, which can be discarded with the filtering pipette. Draw up the filtrate with a Pasteur pipette, and add it to the NMR tube. Add additional CDCl_3 solvent to the NMR tube until the liquid level reaches 50 mm. Determine the ^1H NMR spectrum and interpret the spectrum. This procedure can be applied to other electron-releasing or unreactive compounds such as iodobenzene, 4-iodotoluene, 1-bromo-2-iodobenzene, and 1-bromo-3-iodobenzene. An interesting result is obtained with methyl 2-iodobenzoate in which the methyl ester is converted to the ethyl ester by a transesterification reaction in ethanol during the course of the Sonogashira coupling reaction.


QUESTIONS

1. Draw the structure of the product expected in the following Sonogashira reactions:

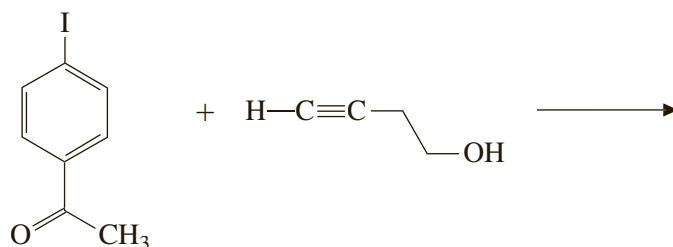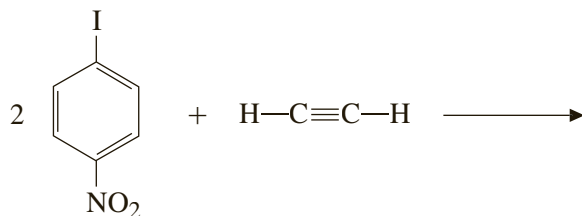

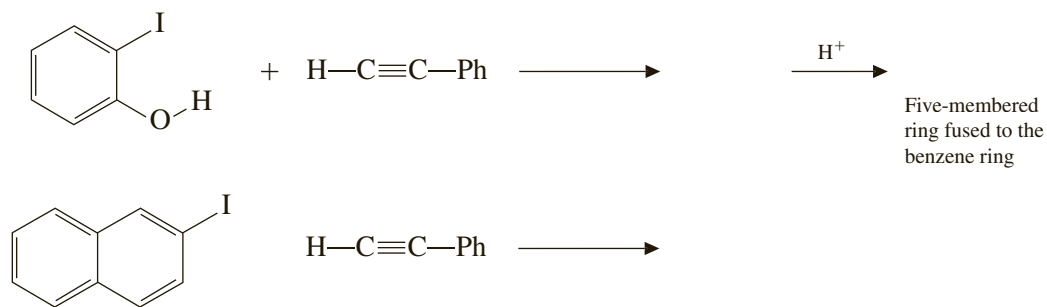

2. Draw the structures of the intermediates and product of the following reaction.

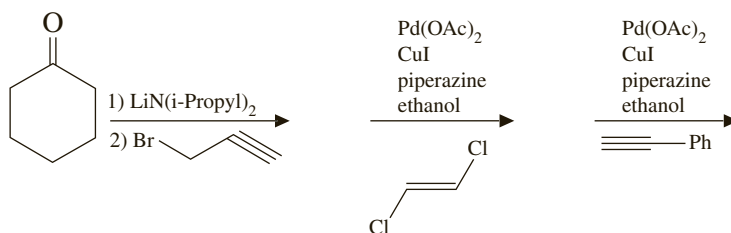

- A small amount of 4,6-decadiyne is formed in reactions involving 1-pentyne. At what point in the mechanism does this compound form?
- Draw a mechanism for the formation of your product.





Grubbs-Catalyzed Metathesis of Eugenol with 1,4-Butenediol to Prepare a Natural Product

Green chemistry


Organometallic chemistry

Ruthenium-catalyzed reactions

Grubbs' catalyst is useful in organometallic chemistry due to its relative stability in air and its tolerance of a variety of solvents. Grubbs' catalyst is a ruthenium-based organometallic catalyst used in cross-coupling metathesis, ring-closing metathesis, and ring-opening metathesis polymerization (ROMP). The four processes are shown below. The dotted line indicates how one can visualize the metathesis process. The development of metathesis reactions in organic synthesis led to the award of the Nobel Prize in Chemistry in 2005 to Yves Chauvin, Robert H. Grubbs, and Richard R. Schrock.

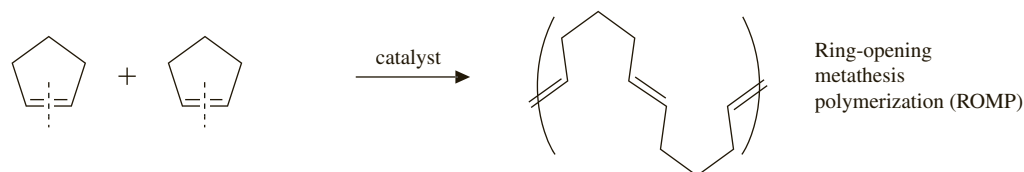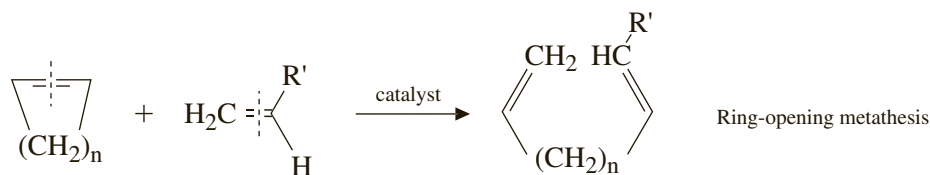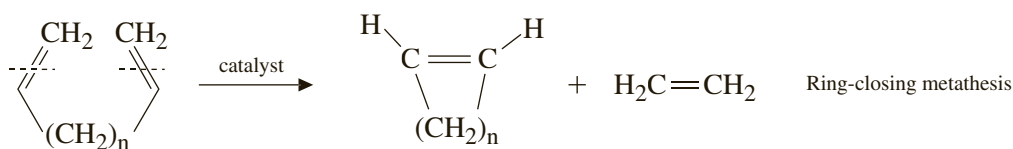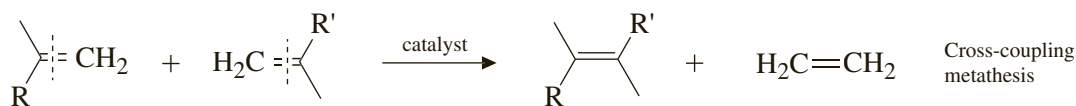

The Grubbs' catalyst that we will use in this experiment is called Grubbs' 2nd Generation catalyst. The IUPAC name is so complicated that researchers don't give the compound a formal IUPAC name! The mechanism for the cross-metathesis reaction is shown on the next page. The current experiment illustrates a very important reaction widely used in research and industry. It is called olefin cross-metathesis.

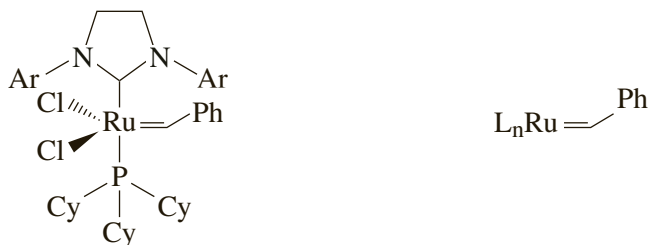

Cy = cyclohexyl

Ph = phenyl


Ar = 2,4,6-trimethylphenyl

Grubbs' Generation 2 catalyst


L_n = ligands

Abbreviated structure for
Grubbs' Generation 2
catalyst

In this experiment, Grubbs' catalyst will be used in the cross-metathesis of eugenol with *cis*-1,4-butenediol¹ to form a natural product known for its medicinal qualities. The product of the reaction, (*E*)-4-(4-hydroxy-3-methoxyphenyl)-2-buten-1-ol, was first isolated from the roots of a South Asian plant, *Zingiber cassumunar*, and is known for its anti-inflammatory and antioxidant properties. You will recognize the pleasant fragrance of eugenol, which is isolated from cloves (see Experiment 15). The reactions are shown below. Natural products such as eugenol are very valuable for making medicinals. The mechanism is shown on the next page.

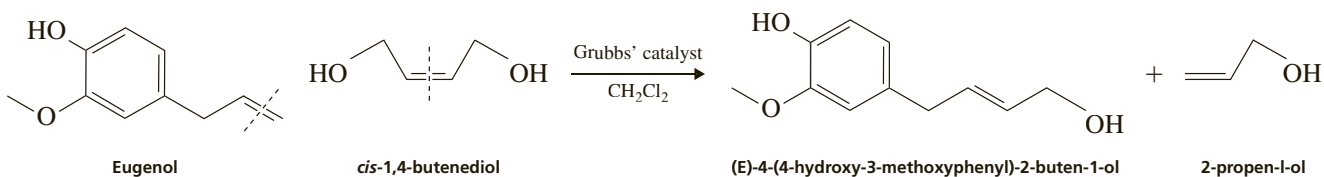

¹ Taber, D. F. and Frankowski, K. J. Grubbs Cross Metathesis of Eugenol with *cis*-1,4-butene-1,4-diol to Make a Natural Product, *Journal of Chemical Education*, **2006**, *83*, 283–284. Experiment developed by Conrardy, D. and Lampman, G. M., Western Washington University, Bellingham, WA.

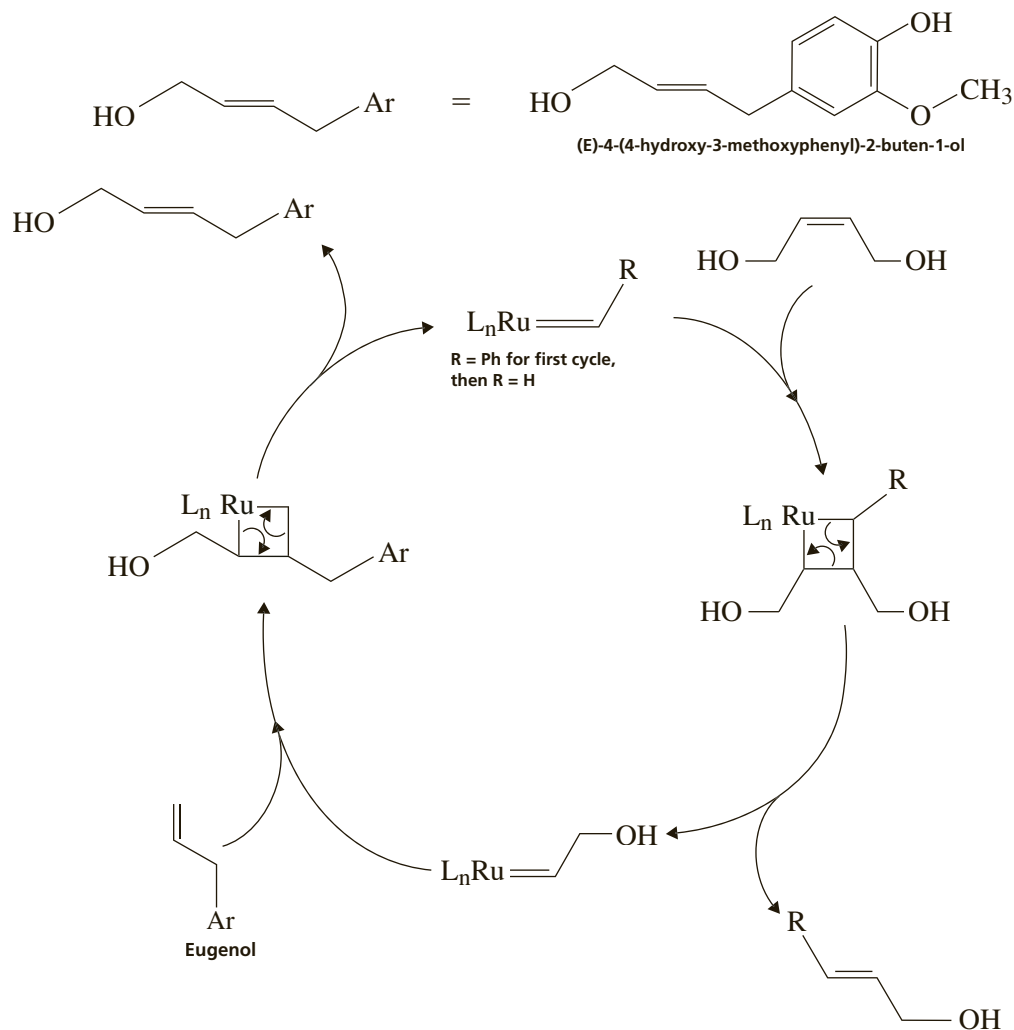

REQUIRED READING

Review: Techniques 5, 6, 7, 12, 19, 26

SPECIAL INSTRUCTIONS

The Grubbs' catalyst is expensive and is air-sensitive. Take care when using it to avoid waste.

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous wastes in the container for aqueous waste. Place the organic waste in the nonhalogenated organic waste container. Place the halogenated waste in the appropriate container.

NOTES TO THE INSTRUCTOR

It is suggested that students work in pairs for this experiment. A rotary evaporator is a useful item to have in the laboratory, but if one is unavailable, the solvents can be removed easily by heating the sample, while blowing air or nitrogen on the sample.

PROCEDURE

Preparation of the Reaction Mixture

Transfer the liquid, eugenol, dropwise to a 50-mL round-bottom flask² using a Pasteur pipette until 0.135 g of eugenol has been obtained. Weigh this material on a 4-place analytical balance. Tare the balance and add 0.490 g of *cis*-1,4-butenediol directly to the same round-bottom flask.

Add 7 mL of methylene chloride to the round-bottom flask and add a small stir bar. Quickly weigh out 0.022 g of Grubbs' 2nd generation catalyst on a piece of weighing paper, using the analytical balance. Weigh the catalyst quickly and add it to the round-bottom flask. The catalyst is sensitive to air and is also very expensive! Work quickly, and remember that it is not important to get an exact amount of the catalyst.

Tightly stopper the flask with a glass stopper or plastic cap to prevent evaporation of the solvent. Cover the cap or glass stopper with Parafilm to reduce the chances of evaporation. Stir the mixture with a magnetic stirrer at a medium rate so as to avoid splashing. If you are using a stirrer/hot plate unit, make sure that the heat is turned off. This reaction proceeds at room temperature. Stir the mixture for at least 1 hour. Allow the mixture to stand at room temperature in your locker, with the stopper or cap securely attached, until the next laboratory period. Allow at least 24 hours. Longer reaction times are also acceptable.

Isolation of the Product

Remove the solvent from the reaction mixture by gently blowing air or nitrogen on the sample while gently heating, or use a rotary evaporator, under vacuum. Continue the evaporation process until a thick, brownish liquid is formed in the bottom of the flask. Remove the flask and add about 1 mL of methylene chloride and about 0.2 g of silica gel³. Swirl the flask so as much of the liquid as possible is absorbed in the silica. Reattach the round-bottom flask to a rotary evaporator and evaporate for another minute or two, under vacuum, to ensure that all of the solvent has been removed. Alternatively, you may blow air or nitrogen on the sample to remove the solvent. A free-flowing solid material will result with the product adsorbed in the silica. Pour the dry solid onto a piece of weighing paper and cover the sample with an inverted beaker.

Column Chromatography

Prepare a silica gel column for chromatography using a 10-mL Pyrex disposable cleanup/drying column (Corning #214210 available from Fisher #05-722-13; the column is about 30 cm long and 1 cm in diameter). Push some cotton down into the bottom using your thermometer. Do not force the cotton too firmly into the tip of the column. It must be tight enough to keep the silica gel from leaking out of the bottom of the column, but not too tight to reduce the flow of solvent. Add enough chromatographic-grade silica gel³ to prepare a 15-cm column.

Make a funnel out of a disposable plastic Pasteur pipette in order to add the sample to the top of the chromatography column. To make the funnel, first cut off the top of a 1-mL plastic pipette and also remove most of the tip to make a small funnel (your instructor should demonstrate this). Pour the silica sample containing

² You may need to obtain this flask from your instructor, along with a glass stopper or plastic cap that fits securely into the standard tapered joint.


³ Fisher Chromatographic Silica Gel, 60-200 mesh #S818-1, DAVISIL® Grade 62, type 150A°.


your adsorbed product from the weighing paper into the top of the silica gel column through the funnel. The solid now resides at the top of the chromatography column.

Add, in portions, 10 mL of petroleum ether (30 to 60°C grade) through the column. Be sure to keep a small amount of liquid at the top of the column at all times to avoid the column drying out. Allow the petroleum ether to flow through the column to wet the silica and begin the elution process. Collect the eluent in an Erlenmeyer flask. Once the petroleum ether has passed through the column, slowly add 30-mL portions of methylene chloride to the column. Allow the column to elute by gravity; do not push the liquid through the column under pressure with a rubber bulb. You are *not likely to see* a distinct band moving down the column; rather, due to dispersion, the colored material spreads out in the column, making it hard to observe the movement of the colored product. The material passing through the column has been variously described as a “trail” of pale light green or a light mint green color or, in some cases, as a grayish/yellow material moving down the column. Because of its pale color, it is often hard to see the material moving down the column. Often the colored material will move below a dark band (you do not want the dark band). Continue to collect the eluent in the Erlenmeyer flask until the colored product reaches the tip of the chromatography column. When the colored product begins to elute, switch from the Erlenmeyer flask to a 50-mL round-bottom flask. You may want to start collecting the eluent early because you may not actually see the colored material dripping out of the tip because the color is so indistinct. If necessary, you may require more methylene chloride to remove the colored product. The liquid in the Erlenmeyer flask is mostly colorless starting material (eugenol), which elutes before the product. The desired product should collect in the round-bottom flask. After all of the colored product has eluted from the column, remove the solvent in the 50-mL round-bottom flask on the rotary evaporator, under vacuum, or blow nitrogen or air to remove solvent.

Isolation and Analysis of the Product

When all of the solvent has been removed, a yellowish-brown solid should be left in the flask; this is the crude product. Add 6 mL of hexane and 1 mL of diethyl ether (not petroleum ether) to the flask and swirl to ensure that all of the product has come into contact with the mixture of solvents. You may need to scrape the bottom of the flask with a spatula to remove the product that is stuck to the bottom. Transfer the product to a Hirsch funnel, under vacuum, to isolate the purified solid product. Use hexane to remove the remaining product from the flask. Continue to draw air through the Hirsch funnel until the product is completely dry. Discard the filtrate. The product should be a solid that ranges in color from yellow to brownish, or perhaps gold or even grayish. Obtain the melting point of the product. Typically, you should expect the melting point to range from 91 to 94°C, but report the actual melting point you obtain. Determine the ^1H NMR spectrum in CDCl_3 . For comparison, the NMR spectrum of the product of the reaction, (*E*)-4-(4-hydroxy-3-methoxyphenyl)-2-buten-1-ol, is shown. The full NMR spectrum is drawn in the lower trace with expansions of individual peaks shown as insets above the full spectrum. The peaks have been labelled on the NMR spectrum to correspond to the structure also shown here.

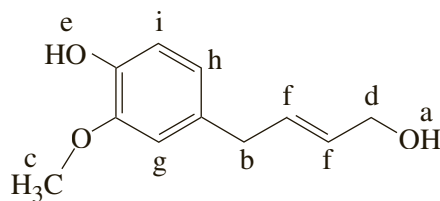


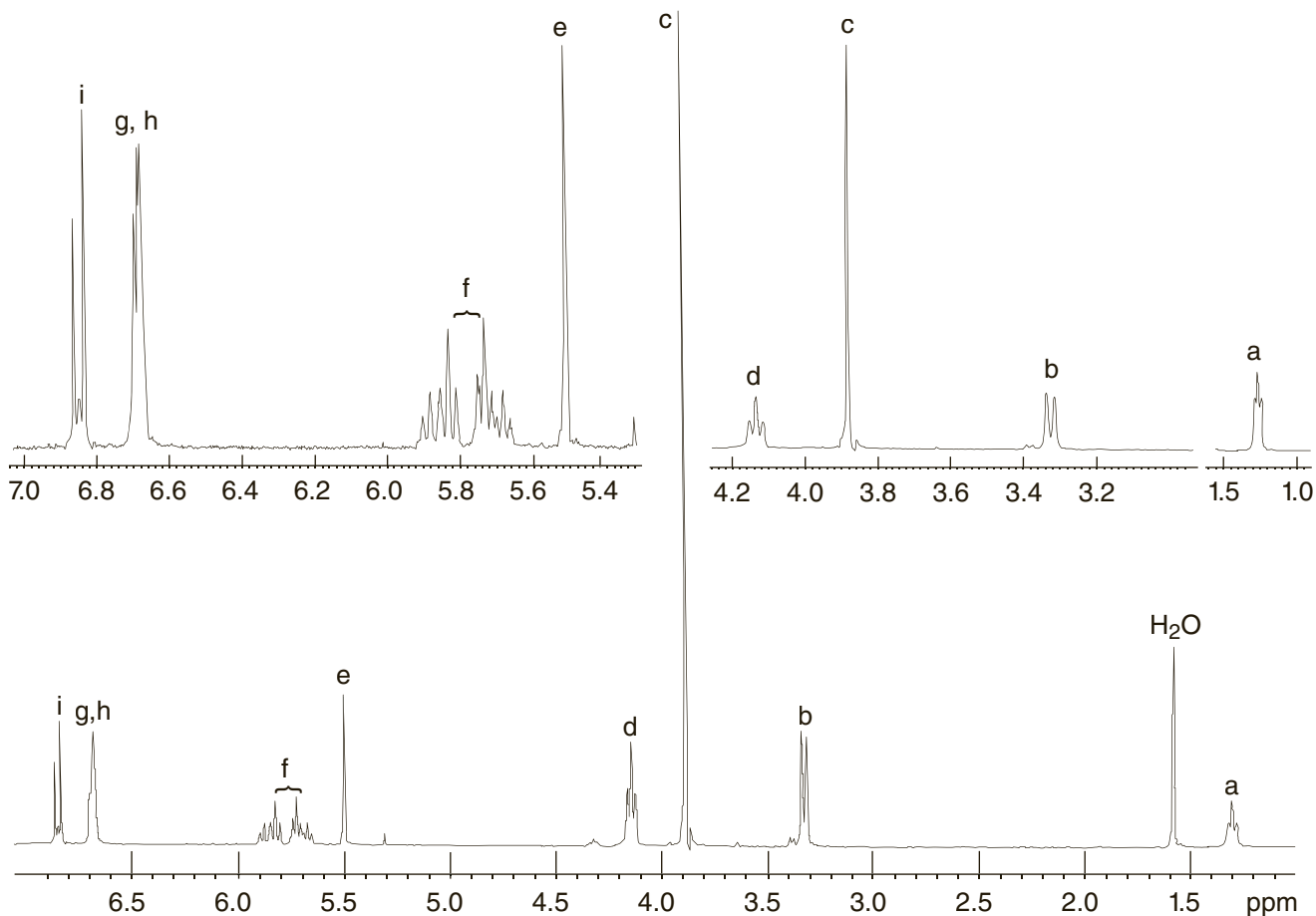



The NMR spectrum of (*E*)-4-(4-hydroxy-3-methoxyphenyl)-2-buten-1-ol. 500 MHz in CDCl₃. The upper traces show expansions for the protons in the metathesis product. The label corresponds to the structure shown above. A peak for water appears at 1.6 ppm.

QUESTIONS

- Column chromatography is used in this experiment to separate the compounds in the mixture from each other. Suggest the order you would expect the following to elute from the column. Use 1 for the first and 4 for the last.
 unreacted eugenol
 unreacted 1,4-butenediol
 ruthenium metal by-products
 your metathesized product
- Draw a mechanism for the following ring-closing metathesis reaction.

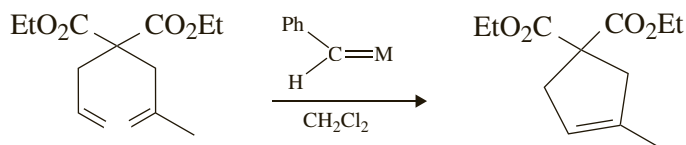

3. Ring-closing metathesis reactions (RCM) have found wide use in forming large ring compounds. Draw the structures of the expected products of the following RCM reactions. See the lab procedure for the general example of RCM.

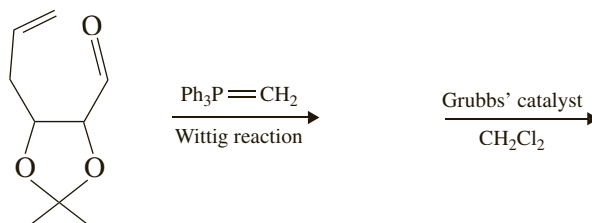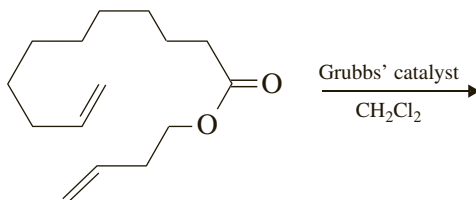

REFERENCES


- Casey, C. P. 2005 Nobel Prize in Chemistry. *J. Chem. Educ.* **2006**, *83*, 192–195.
- France, M. B.; Uffelman, E. S. Ring-Opening Metathesis Polymerization with a Well-Defined Ruthenium Carbene Complex. *J. Chem. Educ.* **1999**, *76*, 661–665.
- Greco, G. E. Nobel Chemistry in the Laboratory: Synthesis of a Ruthenium Catalyst for Ring-Closing Olefin Metathesis. *J. Chem. Educ.* **2007**, *84*, 1995–1997.
- Pappenfus, T. M.; Hermanson, D. L.; Ekerholm, D. P.; Lilliquist, S. L.; Mekoli, M. L. Synthesis and Catalytic Activity of Ruthenium-Indenylidene Complexes for Olefin Metathesis. *J. Chem. Educ.* **2007**, *84*, 1998–2000.
- Scheiper, B.; Glorius, F.; Leitner, A.; Fürstner, A. Catalysis-Based Enantioselective Total Synthesis of the Macrocyclic Spermidine Alkaloid Isoconcinotin. *Proc. Natl. Acad. Sci.* **2004**, *101*, 11960–11965.
- Scholl, M.; Ding, S.; Lee, C. W.; Grubbs, R. H. Synthesis and Activity of a New Generation of Ruthenium-Based Olefin Metathesis Catalysts Coordinated with 1,3-dimethyl-4,5-dihydroimidazol-2-ylidene Ligands. *Org. Lett.* **1999**, 953–956.
- Schrock, R. R. Living Ring-Opening Metathesis Polymerization Catalyzed by Well-Characterized Transition-Metal Alkylidene Complexes. *Acc. Chem. Res.* **1990**, *23*, 158–165.
- Taber, D. F.; Frankowski, K. J. Grubbs Cross Metathesis of Eugenol with *cis*-1,4-Butene-1, 4-diol to Make a Natural Product. *J. Chem. Educ.* **2006**, *83*, 283–284.
- Trnka, T. M.; Grubbs, R. H. Development of L2X2Ru=CHR Olefin Metathesis Catalysts: An Organometallic Success Story. *Acc. Chem. Res.* **2001**, *34*, 18–29.

The Aldol Condensation Reaction: Preparation of Benzalacetophenones (Chalcones)

Aldol condensation

Crystallization

Benzaldehyde reacts with a ketone in the presence of base to give α, β -unsaturated ketones. This reaction is an example of a crossed aldol condensation where the intermediate undergoes dehydration to produce the resonance-stabilized unsaturated ketone.


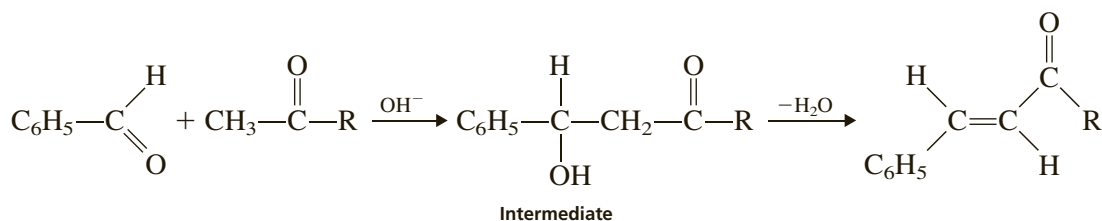


Crossed aldol condensations of this type proceed in high yield because benzaldehyde cannot react with itself by an aldol condensation reaction because it has no α -hydrogen. Likewise, ketones do not react easily with themselves in aqueous base. Therefore, the only possibility is for a ketone to react with benzaldehyde.

In this experiment, procedures are given for preparing benzalacetophenones (chalcones). You should choose one of the substituted benzaldehydes and react it with the ketone, acetophenone. All the products are solids that can be recrystallized easily.

Benzalacetophenones (chalcones) are prepared by the reaction of a substituted benzaldehyde with acetophenone in aqueous base. Piperonaldehyde, *p*-anisaldehyde, and 3-nitrobenzaldehyde are used.

An optional molecular modeling exercise is provided in this experiment. We will examine the reactivity of the enolate ion of a ketone to see which atom, oxygen or carbon, is more nucleophilic. The molecular modeling part of this experiment will help you to rationalize the experimental results of this experiment. It would be helpful to look at Experiment 20E, in addition to the material given in this experiment.

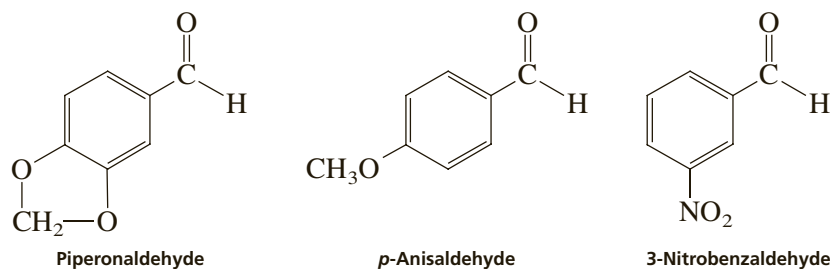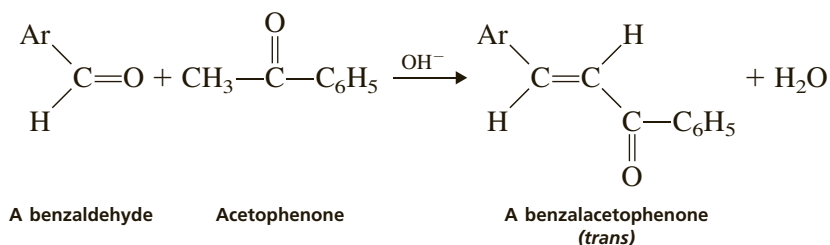

REQUIRED READING

Review: Technique 8 Section 8.7
Technique 11 Sections 11.3 and 11.4

SPECIAL INSTRUCTIONS

Before beginning this experiment, select one of the substituted benzaldehydes. Alternatively, your instructor may assign a particular compound to you.

SUGGESTED WASTE DISPOSAL

All filtrates should be poured into a waste container designated for nonhalogenated organic waste.

PROCEDURE

Running the Reaction

Choose one of three aldehydes for this experiment: piperonaldehyde (solid), 3-nitrobenzaldehyde (solid), or *p*-anisaldehyde (liquid). Place 0.150 g of piperonaldehyde (3,4-methylenedioxybenzaldehyde, *MW* = 150.1) or 0.151 g of 3-nitrobenzaldehyde (*MW* = 151.1) into a 5-mL conical vial. Alternatively, transfer 0.13 mL of *p*-anisaldehyde (4-methoxybenzaldehyde, *MW* = 136.2) to a *tared* conical vial and reweigh the vial to determine the weight of material transferred.

Add 0.12 mL of acetophenone (*MW* = 120.2, *d* = 1.03 g/mL) and 0.80 mL of 95% ethanol to the vial containing your choice of aldehyde. Place the conical vial into a 50-mL beaker. Stir the mixture with a microspatula to dissolve any solids present. You may need to warm the mixture on a hot plate to dissolve the solids. If this is necessary, then cool the solution to room temperature before proceeding with the next step.

Add 0.10 mL of aqueous sodium hydroxide solution¹ to the aldehyde/ acetophenone mixture. Stir the mixture with your microspatula until it solidifies. Before the mixture solidifies, you may observe some cloudiness. *Wait until the cloudiness has been replaced with an obvious precipitate settling out to the bottom of the conical vial before proceeding to the next paragraph.* Continue stirring with your microspatula until a solid forms (approximately 3 to 5 minutes).² Scratching the inside of the conical vial with your microspatula may help to crystallize the chalcone.

Isolation of the Crude Chalcone

Add 2 mL of ice water to the vial *after a solid has formed, as indicated in the previous paragraph.* Stir the solid in the mixture with a spatula to break up the solid mass.

Transfer the mixture to a small beaker with 3 mL of ice water. Stir the precipitate to break it up and then collect the solid on a Hirsch funnel. Wash the product with cold water. Let the solid air-dry for about 30 minutes. Weigh the solid and determine the percentage yield.

Crystallization of the Chalcone

Purify all of the crude chalcone from hot 95% ethanol or hot methanol using the semimicroscale crystallization procedure (Technique 11, Section 11.3). Alternatively, purify part of the chalcone using a Craig tube (Technique 8, Section 8.7, and Technique 11, Section 11.4), as follows:

3,4-Methylenedioxychalcone (from piperonaldehyde) Crystallize a 0.040-g sample from about 0.5 mL of hot 95% ethanol; literature melting point is 122°C.

4-Methoxychalcone (from p-anisaldehyde) Crystallize a 0.075-g sample from about 0.3 mL of hot 95% ethanol. Scratch the tube to induce crystallization while cooling; literature melting point is 74°C.

3-Nitrochalcone (from 3-nitrobenzaldehyde) Crystallize a 0.025-g sample from about 1 mL of hot methanol. Scratch the tube gently to induce crystallization while cooling; literature melting point is 146°C.

Laboratory Report


Determine the melting point of your purified product. At the option of the instructor, obtain the proton and/or ¹³C NMR spectrum. Include a balanced equation for the reaction in your report. Submit the crude and purified samples to the instructor in labeled vials.

MOLECULAR MODELING (OPTIONAL)

In this exercise we will examine the enolate ion of acetone and determine which atom, oxygen or carbon, is the more nucleophilic site. Two resonance structures can be drawn for the enolate ion of acetone, one with the negative charge on oxygen, structure **A**, and one with the negative charge on carbon, structure **B**.

¹ The instructor should prepare the concentrated aqueous sodium hydroxide in advance, in the ratio of 0.60 g of sodium hydroxide to 1 mL water.

² In some cases, the chalcone may not precipitate. If this is the case, cap the conical vial and allow it to stand until the next laboratory period. Usually, the chalcone will precipitate during that time. An additional portion of base will sometimes be helpful, as will *gentle* warming.

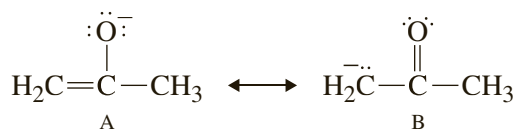



The enolate ion is an **ambident nucleophile**, a nucleophile that has two possible nucleophilic sites. Resonance theory indicates that structure **A** should be the major contributing structure because the negative charge is better accommodated by oxygen, a more electronegative atom than carbon. However, the reactive site of this ion is carbon, not oxygen. Aldol condensations, brominations, and alkylations take place at carbon, not oxygen. In frontier molecular orbital terms (see the essay “Computational Chemistry-ab Initio and Semiempirical Methods”), the enolate ion is an electron pair donor, and we would expect the pair of electrons donated to be those in the highest occupied molecular orbital, the HOMO.





In the structure-building editor of your modeling program, build structure **A**. Be sure to delete an unfilled valence from oxygen and to place a -1 charge on the molecule. Request a geometry optimization at the AM1 semiempirical level. Also request the HOMO surface and maps of the HOMO and the electrostatic potential onto the electron density surface. Submit your selections for computation. Plot the HOMO on the screen. Where are the biggest lobes of the HOMO, on carbon or on oxygen? Now map the HOMO onto the electron density surface. The “hot spot,” the place where the HOMO has the highest density at the point where it intersects the surface, will be bright blue. What do you conclude from this mapping? Finally, map the electrostatic potential onto the electron density. This shows the electron distribution in the molecule. Where is the overall electron density highest, on oxygen or on carbon?

Finally, build structure **B** and calculate the same surfaces as requested for structure **A**. Do you obtain the same surfaces as for structure **A**, or are they different? What do you conclude? Include your results, along with your conclusions, in your report on this experiment.


QUESTIONS

1. Give a mechanism for the preparation of the appropriate benzalacetophenone using the aldehyde and ketone that you selected in this experiment.
2. Draw the structure of the *cis* and *trans* isomers of the compound that you prepared. Why did you obtain the *trans* isomer?
3. Using infrared spectroscopy, how could you experimentally determine that you have the *trans* isomer rather than the *cis* one? (*Hint*: See Technique 25, Section 25.14.)
4. Provide the starting materials needed to prepare the following compounds:

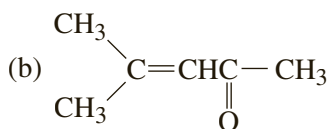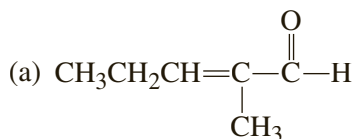

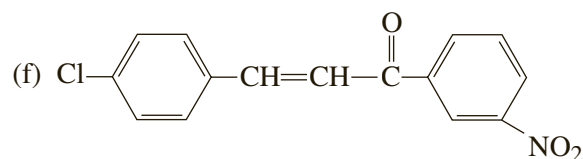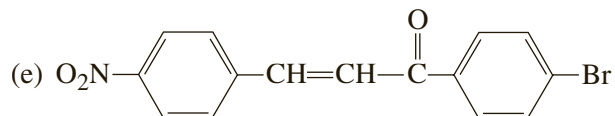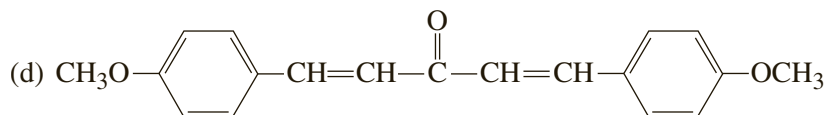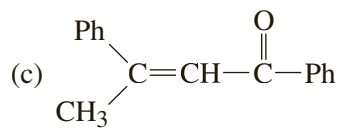

5. Prepare the following compounds starting from benzaldehyde and the appropriate ketone. Provide reactions for preparing the ketones starting from aromatic hydrocarbon compounds (see Experiment 60).

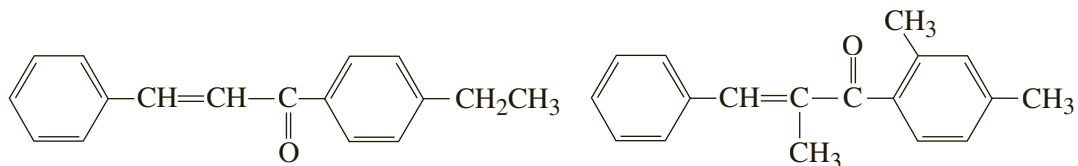


A Green Enantioselective Aldol Condensation Reaction

Green chemistry

Proline-catalyzed asymmetric induction


The aldol condensation is a fundamental reaction in chemistry and biology. In its most common form, a ketone reacts with an aldehyde to form a 3-hydroxy ketone (sometimes referred to as a β -hydroxy ketone). A new C—C bond is formed in the reaction, and a new stereocenter is formed at the position of the hydroxyl group. The most common catalyst used in aldol condensation reactions is sodium hydroxide. Under these conditions, a racemic mixture is formed when acetone is allowed to react with an aldehyde. In the example shown, acetone is reacted with isobutyraldehyde.

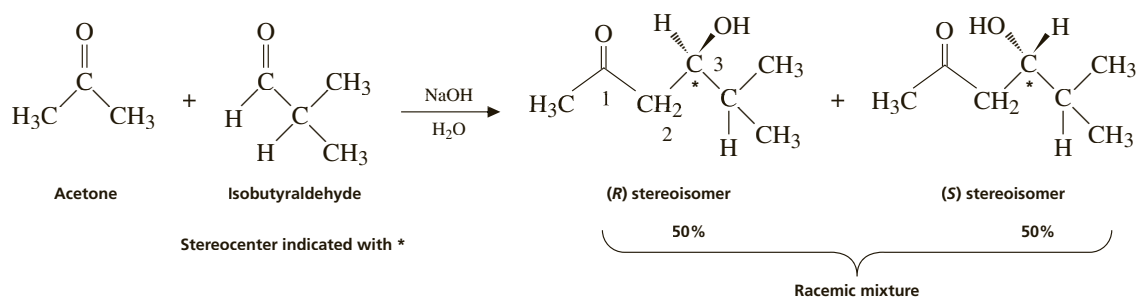

The dream of every organic chemist is to avoid creating a racemic mixture and instead obtain a single stereoisomer! This type of reaction is often referred to as an enantioselective reaction, in which one stereoisomer is primarily created in the reaction. In order to do this, one needs to start with a chiral catalyst, in this case *L*-proline, a naturally occurring amino acid. Biological reactions form one stereoisomer because the enzymes in natural systems are themselves chiral. In effect, we are trying to mimic the process that occurs in natural systems. *L*-proline mimics the class I aldolase enzymes in natural systems.^{1, 2}


¹Bennett, G. D. A Green Enantioselective Aldol Condensation for the Undergraduate Organic Laboratory. *J. Chem. Educ.* **2006**, *83*, 1871–1872. Experiment developed by Bowen, G. and Lampman, G. M., Western Washington University, Bellingham WA.

²List, B.; Lerner, R. A.; Barbas III, C. F. Proline-Catalyzed Direct Asymmetric Aldol Reactions. *J. Am. Chem. Soc.* **2000**, *122*, 2395–2396.

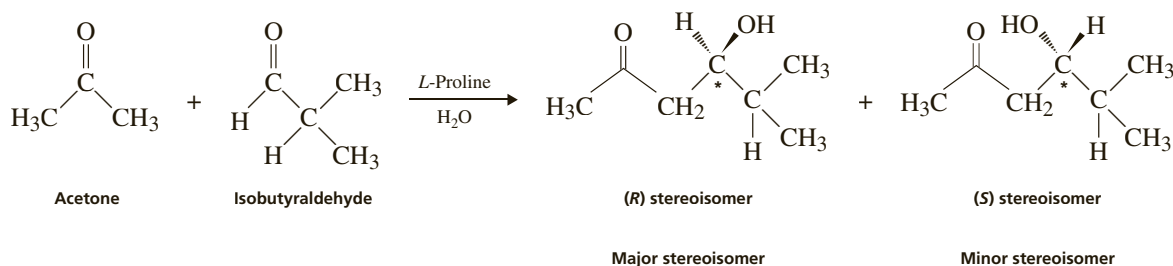

This experiment demonstrates an important concept that has wide use in the pharmaceutical industry, where formation of single enantiomers is critical. In many cases, one enantiomer elicits the correct biological response, while the other enantiomer may have harmful effects.

Often the product from the aldol condensation reaction undergoes further reaction by eliminating the elements of water. This is especially common when a substituted acetophenone reacts with substituted benzaldehyde. Two experiments are included in this book to demonstrate this pathway (see Experiment 39 and 65). In these types of experiments, the intermediate β -hydroxyketone loses water to form a conjugated ketone. The driving force for this reaction is the formation of the highly resonance-stabilized ketone. The compounds formed in this reaction are given the trivial name of chalcone or the IUPAC name of 1,3-diphenyl-2-propen-1-one. The chalcone that is formed loses the stereocenters and becomes achiral (non-chiral). We should expect that in certain reactions, the aldol product will give rise to some elimination by-product. Fortunately, elimination is not a major pathway for the reaction of acetone with isobutyraldehyde.

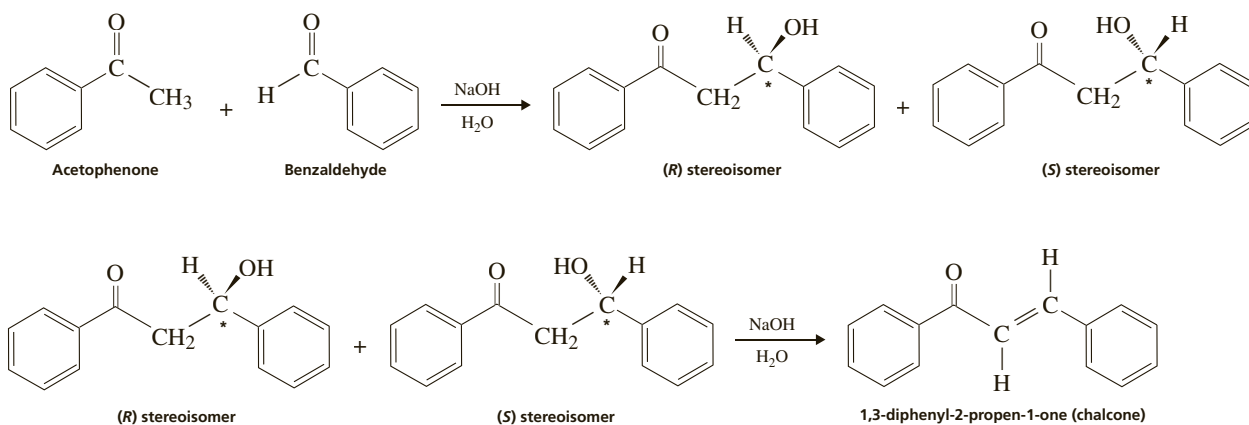

The relative amounts of the aldol condensation product and the elimination (dehydration) products can be determined by NMR. We will employ a polarimeter to determine the degree of stereospecificity for the reaction of acetone with isobutyraldehyde to give the aldol adduct. In order to determine the enantiomeric excess (ee), we need the specific rotation value for one of the pure enantiomers, in this case the (R) enantiomer. Unfortunately, it is sometimes difficult to find the specific

rotation values for other reactions in the chemical literature.³ Other methods must be employed in research laboratories to determine the enantioselectivity of the aldol condensation products. Methods that are often used in research are chiral gas chromatography and chiral HPLC.

The mechanism for the *L*-proline-catalyzed aldol condensation reaction is shown in Scheme 1. Steps 3A and 3B show the two possible chair structures for the transition states leading to the (*R*) and (*S*) products. Notice that the isopropyl group is in an axial position in 3A, while this group is attached to an equatorial position in 3B. We should, therefore, expect that the lower energy barrier should proceed through 3B and yield the (*R*) adduct as the major aldol condensation product.

REQUIRED READING

Review: Techniques 5, 6, 7, 8, 12, 25, and 26

New: Technique 23 Polarimetry

SPECIAL INSTRUCTIONS

Isobutyraldehyde is an irritant, and some of the products can cause an allergic response. It is advised that you wear gloves for this experiment.

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous wastes in the container for aqueous waste. Place the organic waste in the nonhalogenated organic waste container.

NOTES TO THE INSTRUCTOR

Although this experiment is labeled as a green experiment, acetone is used in a large excess—making the experiment rather poor in terms of atom economy, even though the addition reaction itself has a high degree of atom economy. The excess is required to avoid unfavorable side reactions. The use of *L*-proline, a natural amino acid, in catalytic amounts helps to make this experiment green. Since the reaction proceeds slowly, it would be difficult to apply this reaction in an industrial setting. Diethyl ether is used as a solvent for extraction, avoiding the use of methylene chloride.

Other substrates may be used in this experiment. Examples include the reaction of acetone with pivaldehyde (2,2-dimethylpropanal) and the reaction of acetone with acetophenone. In the latter case, a significant amount of elimination is observed. The NMR easily reveals when elimination is an important side reaction. NMR itself cannot be used to determine the relative amounts of the enantiomers, but it is the best way of analyzing the amount of elimination that occurs in these

³ Ramachandran, P. V.; Xu, Wei-chu; Brown, H. C. Contrasting Steric Effects of the Ketones and Aldehydes in the Reactions of the Diisopinocampheyl Enolborinates of Methyl Ketones with Aldehydes. *Tetrahed. Lett.* **1996**, 37, 4911–4914.

reactions. For example, you can expect about 3% to 6% elimination (dehydration) in the acetone–isobutyraldehyde or acetone–pivaldehyde reactions. The acetone–benzaldehyde reaction gives more than 60% of the elimination product.

Unfortunately, when the adducts are analyzed by GC-MS, even more of the elimination (dehydration) product forms in the heated inlet of the gas chromatograph. Therefore, it is recommended that the NMR be used to determine the relative amounts of the aldol adduct and dehydration products. The polarimeter is used to determine the ee in the *L*-proline-catalyzed reaction.

Using polarimetry, you may find that the class will obtain a value of $+34^\circ$ for the adduct formed by acetone and isobutyraldehyde. This value is compared to the specific rotation for the pure (*S*) enantiomers of 61.7° to give an enantiomeric excess value of 55%. It has been suggested that lower values occur when traces of water are present during the course of the *L*-proline-catalyzed reaction. The acetone–pivaldehyde reaction gives an adduct that has a specific rotation of 56.9° , which yields a calculated value of 69% for the enantiomeric excess.⁴

PROCEDURE

Transfer 1.0 mL of isobutyraldehyde to a preweighed 25-mL round-bottom flask using a 1000 μ L automatic pipette, and reweigh the flask to determine the precise weight of isobutyraldehyde transferred to the flask. Add 14 mL of acetone and 0.23 g (2 mmol) of *L*-proline to the flask. Add a magnetic stir bar and insert a glass stopper or plastic cap into the neck of the flask. Stir the mixture for 1 week at room temperature.

Pour the contents of the round-bottom flask into a beaker. Add 20 mL of diethyl ether to the beaker. Add another 5 mL of diethyl ether to rinse out the round-bottomed flask. Some solid may be present that does not dissolve (discard it). Pour 50 mL of saturated aqueous sodium chloride solution into the beaker. Transfer all of the diethyl ether and the saturated salt solution into a separatory funnel, avoiding adding the stir bar to the separatory funnel. Shake the funnel in order to ensure that the product is extracted into the diethyl ether. Drain the lower aqueous layer and discard it. Pour the diethyl ether extract from the top of the separatory funnel into an Erlenmeyer flask and add anhydrous magnesium sulfate to dry the extract. Remove the drying agent by gravity filtration through a fluted filter into a preweighed round-bottom flask. Remove the solvent with a rotary evaporator, under vacuum. Attach the flask to a good vacuum pump to remove any remaining acetone and diethyl ether. Weigh the flask to determine the yield of aldol condensation product and calculate the percentage yield for the reaction.

Add the product to a preweighed 5-mL volumetric flask, and reweigh the flask after the addition. Dissolve the sample in chloroform up to the mark on the volumetric flask. Calculate the density in g/mL for the chloroform solution. Add the chloroform solution to a 0.5-dm polarimeter cell and obtain the optical rotation for the sample. Calculate the specific rotation using the equation shown in Technique 23, Section 23.2. The optically pure (*R*) enantiomer has a reported specific rotation

⁴Ramachandran, P. V.; Xu, Wei-chu; Brown, H. C. Contrasting Steric Effects of the Ketones and Aldehydes in the Reactions of the Diisopinocampheyl Enolborinates of Methyl Ketones with Aldehydes. *Tetrahed. Lett.* **1996**, 37, 4911–4914.


value of $+61.7^\circ$.⁵ The ee is obtained using the equation shown in Technique 23, Section 23.5.⁶ Use the ee value to calculate the percentages of the (*S*) and (*R*) enantiomers in the mixture (see Technique 23, Section 23.5).

At the option of your instructor, determine the infrared spectrum and the high field (500 MHz) ^1H NMR spectrum for your aldol condensation product, (*R*)-4-hydroxy-5-methyl-2-hexanone. Figure 1 shows the ^1H NMR spectrum of the product. All of the peaks are assigned to the structure shown, except for the OH group. The stereocenter in the product introduces some very interesting features to the NMR spectrum shown in Figure 1! Of particular interest in the NMR spectrum is the area between 2.50 to 2.65 ppm in the product. This area reveals the presence of the two nonequivalent diastereotopic protons in the methylene group, H_a and H_b , shown in expansion in Figure 1. The peaks in this expansion are labeled with Hz values so that coupling constants can be calculated.

H_a centers on 2.62 ppm and is a doublet of doublets, yielding a coupling constant, $^2J_{ab} = 17.5$ Hz (1317.38 – 1299.80 Hz). In addition, H_a is coupled to H_c , yielding a value for $^3J_{ac} = 2.4$ Hz (1302.25 – 1299.80 Hz). The other diastereotopic proton, H_b , centers on about 2.54 ppm and is also a doublet of doublets, with $^2J_{ab} = 17.5$ Hz (1281.25 – 1263.67 Hz) and $^3J_{bc} = 9.7$ Hz (1281.25 – 1271.48 Hz). Notice that the diastereotopic protons, H_a and H_b , have identical 2J values of 17.5 Hz. The 3J values are different because the dihedral angles are not the same. The dihedral angle for protons H_a and $\text{H}_c = 60^\circ$, whereas the angle for protons H_b and $\text{H}_c = 180^\circ$. To summarize, $^2J_{ab} = 17.5$ Hz, $^3J_{bc} = 9.7$ Hz, and $^3J_{ac} = 2.4$ Hz.

The other area of interest in the ^1H NMR spectrum is the pair of doublets, labeled as f on the structure and spectrum, that appear at 0.916 and 0.942 ppm in the expansion in Figure 1. It turns out that the two methyl groups are also nonequivalent because of the presence of the stereocenter, and are also diastereotopic. Because of the presence of the stereocenter, the methyl groups appear at different places in the NMR spectrum.

We cannot determine the percentages of the two possible enantiomers in the NMR spectrum shown here. It should be strongly stated that both the (*R*) and the (*S*) enantiomers will have identical NMR spectra! *Only if the two enantiomers are placed in a chiral environment will they have different NMR spectra!* A polarimeter will show different behavior for the two enantiomers because there is a chiral environment present.

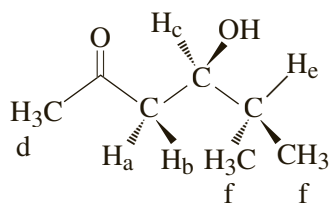

⁵ List, B.; Lerner, R. A.; Barbas III, C. F. Proline-Catalyzed Direct Asymmetric Aldol Reactions. *J. Am. Chem. Soc.* **2000**, *122*, 2395–2396. These researchers reported an ee as high as 96% for the reaction of isobutyraldehyde and acetone with *L*-proline.

⁶ A typical value of 60% ee may be obtained yielding a mixture of 80% (*R*) and 20% (*S*) enantiomers. Other methods were attempted in an effort to obtain the ee: chiral GC column and chiral chemical shift reagents (see Technique 26, Section 26.15), without success. Better results may be achieved by using more concentrated samples and a smaller cell so that higher values of the rotation may be obtained.

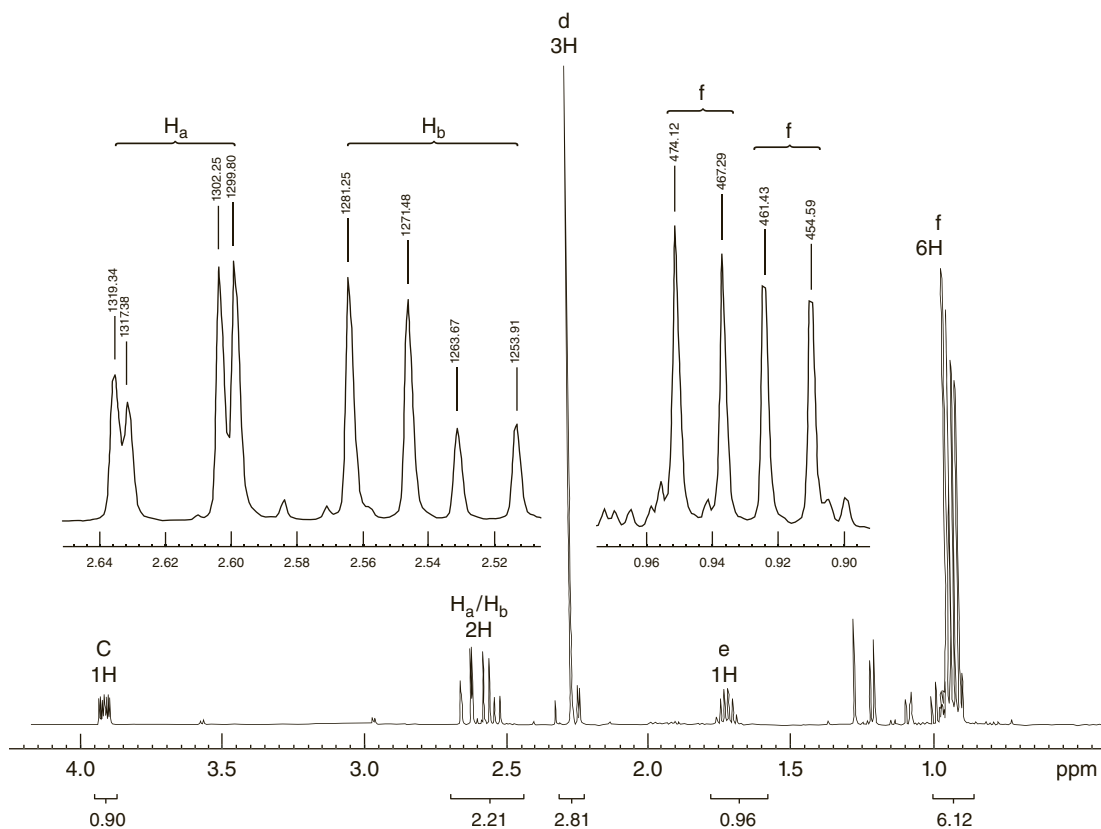
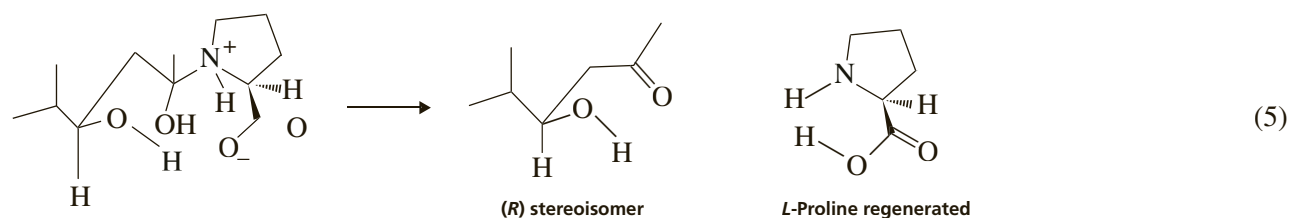
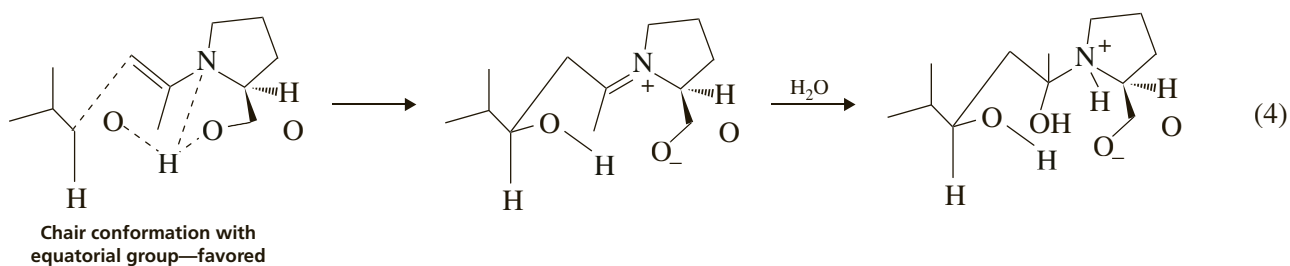
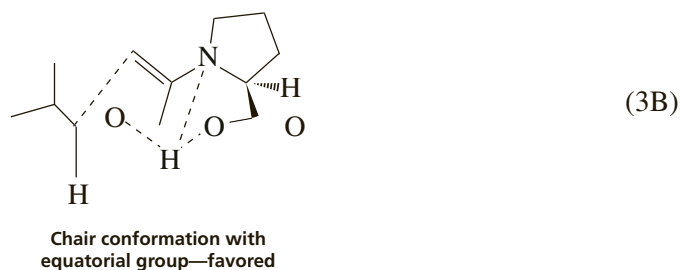
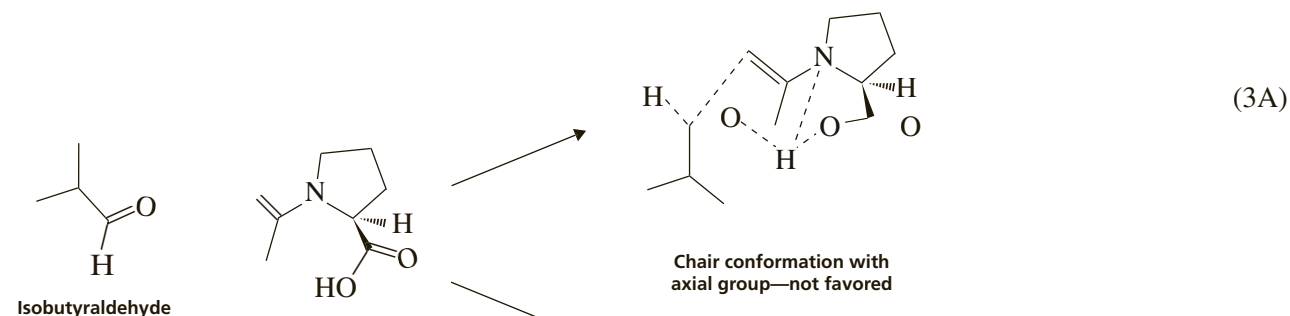
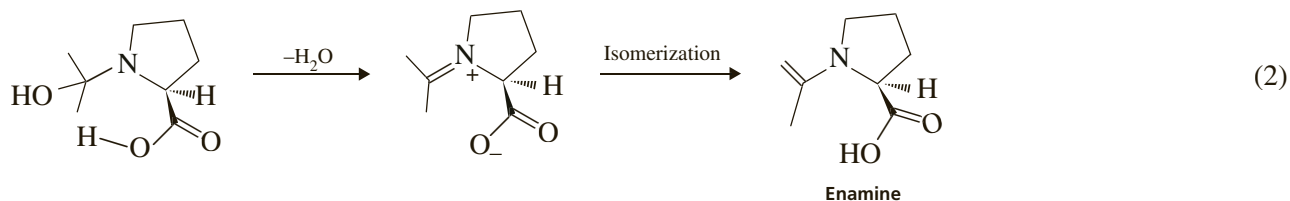
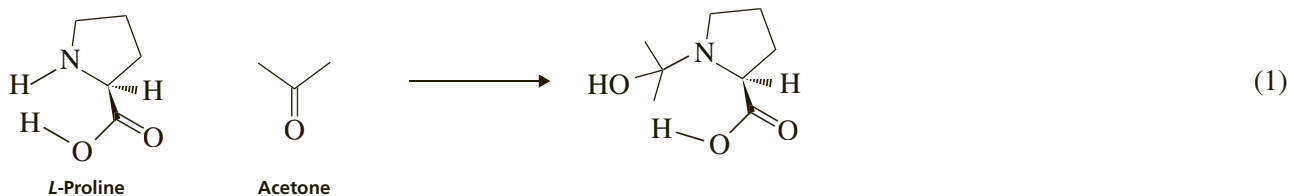


Figure 1. 500 MHz ^1H NMR spectrum of the *L*-proline-catalyzed aldol condensation of isobutyraldehyde and acetone. The upper traces show expansions of the diastereotopic methylene group, H_a and H_b , appearing between 2.50 to 2.65 ppm. Also shown in the upper traces are expansions for the two diastereotopic methyl groups, labeled as f on the spectrum. Impurity peaks appear between 0.96 to 1.3 ppm.

Scheme 1. Mechanism of the *L*-proline catalyzed aldol condensation of isobutyraldehyde and acetone.


Preparation of an α,β -Unsaturated Ketone via Michael and Aldol Condensation Reactions

Crystallization


Michael reaction (conjugate addition)

Aldol condensation reaction

This experiment illustrates how two important synthetic reactions can be combined to prepare an α,β -unsaturated ketone, 6-ethoxycarbonyl-3,5-diphenyl-2-cyclohexenone. The first step in this synthesis is a sodium hydroxide-catalyzed conjugate addition of ethyl acetoacetate to *trans*-chalcone (a Michael addition reaction). Sodium hydroxide serves as a source of hydroxide ion to catalyze the reaction.¹ In the reactions that follow, Ph and Et are abbreviations for the phenyl and ethyl groups, respectively.

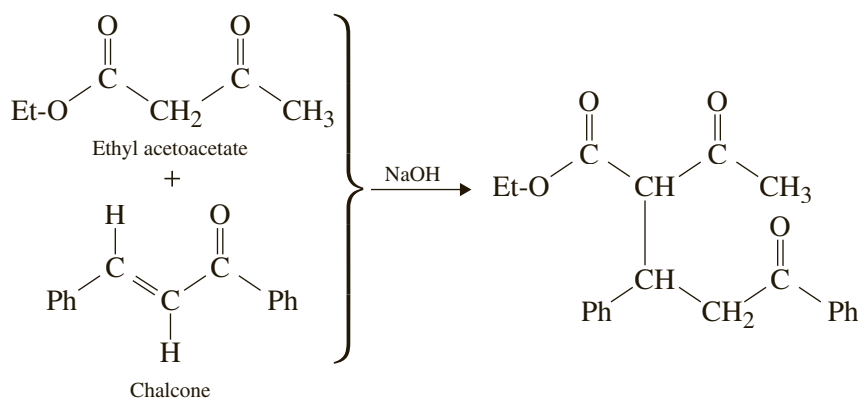


The second step of the synthesis is a base-catalyzed aldol condensation reaction. The methyl group loses a proton in the presence of base, and the resulting methylene carbanion nucleophilically attacks the carbonyl group. A stable six-membered ring is formed. Ethanol supplies a proton to yield the aldol intermediate.

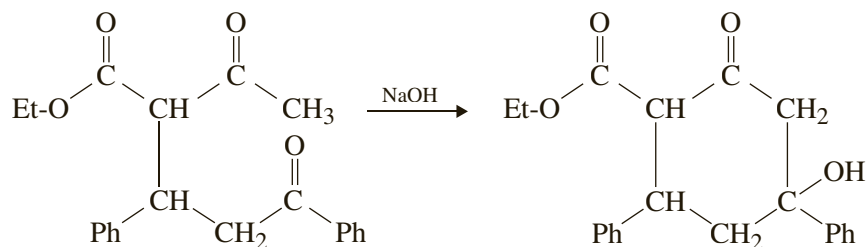

¹ Barium hydroxide has also been used as a catalyst (see References).

Finally, the aldol intermediate is dehydrated to form the final product, 6-ethoxycarbonyl-3,5-diphenyl-2-cyclohexenone. The α,β -unsaturated ketone that is formed is stable because of the conjugation of the double bond with both the carbonyl group and a phenyl group.

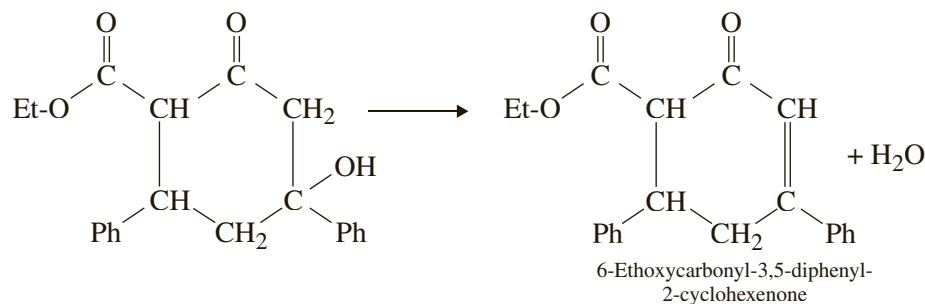

REQUIRED READING

Review: Techniques 7, 8, 11, 12

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous wastes containing ethanol in the bottle designated for aqueous wastes. Ethanolic filtrates from the crystallization of the product should be poured into the nonhalogenated organic waste container.

NOTES TO THE INSTRUCTOR

The *trans*-chalcone (Aldrich Chemical Co., #13,612-3) should be finely ground for use by the class.

PROCEDURE

Assembling the Apparatus

To a 10-mL round-bottom flask, add 0.24 g of finely ground *trans*-chalcone, 0.15 g of ethyl acetoacetate, and 5 mL of absolute ethanol. Swirl the flask until all or most of the solid dissolves and place a boiling stone in the flask. Add 0.25 mL of 2.2M NaOH to the mixture. Attach a water-jacketed condenser to the round-bottom flask and heat the mixture to reflux using an aluminum block and hot plate. Once the mixture has been brought to a gentle boil, continue to reflux the mixture for at least 1 hour. During this reflux, the mixture will become very cloudy.

Isolation of the Crude Product

After the end of the reflux period, let the mixture cool to room temperature. Add 2 mL of water and scratch the inside of the flask with a glass stirring rod to induce crystallization (an oil may form; scratch vigorously). Place the flask in an ice bath for a minimum of 30 minutes. It is essential to cool the mixture thoroughly in order to completely crystallize the product. Because the product may precipitate slowly,

you should also scratch the inside of the flask occasionally over the 30-minute period, as well as cool it in an ice bath.

Vacuum-filter the crystals on a Hirsch funnel, using 1 mL of ice-cold water to aid in the transfer. Then rinse the round-bottom flask with 1 mL of ice-cold 95% ethanol to complete the transfer of the remaining solid from the flask to the Hirsch funnel. Allow the crystals to air-dry overnight. Alternatively, the crystals may be dried for 30 minutes in an oven set at 75–80°C. Weigh the dry product. The solid contains some sodium hydroxide and sodium carbonate, which are removed in the next step.

Removal of Catalyst

Place the solid product in a test tube. Add 1.5 mL of reagent-grade acetone and stir the mixture with a spatula.² Most of the solid dissolves in acetone, but do not expect all of it to dissolve. Using a Pasteur pipette, remove the liquid and transfer it to a glass centrifuge tube, leaving as much solid as possible behind in the test tube. It is impossible to avoid drawing some solid up into the pipette, so the transferred liquid will contain suspended solids and the solution will be very cloudy. You should not be concerned about the suspended solids in the cloudy acetone extract, because the centrifugation step will clear the liquid completely. Centrifuge the acetone extract for approximately 2–3 minutes or until the liquid clears. Using a clean, dry Pasteur pipette, transfer the *clear* acetone extract from the centrifuge tube to a *dry, preweighed* test tube. If the transfer operation is done carefully, you should be able to leave the solid behind in the centrifuge tube. The solids left behind in the test tube and centrifuge tube are inorganic materials related to the sodium hydroxide originally used as the catalyst.

Evaporate the acetone solvent by carefully heating the test tube in a hot water bath while directing a light stream of dry air or nitrogen in the tube. Use a *slow* stream of gas to avoid blowing your product out of the tube. When the acetone has evaporated, you may be left with an oily solid in the bottom of the tube. Scratch the oily product with a spatula to induce crystallization. You may need to redirect air or nitrogen in the test tube to remove all traces of acetone. Reweigh the test tube to determine the yield of this partially purified product.

Crystallization of Product

Crystallize the product in a 10-mL Erlenmeyer flask using a minimum amount (approximately 2 mL) of hot 95% ethanol.³ After the solid has dissolved, allow the flask to cool slightly. Scratch the inside of the flask with a glass stirring rod until crystals appear. Allow the flask to sit undisturbed at room temperature for a few minutes. Then place the flask in an ice-water bath for at least 15 minutes.

Collect the crystals by vacuum filtration on a Hirsch funnel. Use two 0.5-mL portions of ice-cold 95% ethanol to aid in the transfer. Allow the crystals to dry overnight or dry them for 30 minutes in a 75–80°C oven. Weigh the dry 6-ethoxycarbonyl-3,5-diphenyl-2-cyclohexenone and calculate the percentage yield. Determine the melting point of the product (literature value, 111–112°C). Submit the sample to the instructor in a labeled vial.

² You may need to add more acetone than indicated in the procedure because a larger yield of product may have been obtained. You can add up to 3 mL of acetone. Excess acetone will not affect the results.

³ The 2 mL of ethanol indicated in the procedure is an approximation. You may need to add *more hot or less hot* 95% ethanol to dissolve the solid. Add boiling ethanol until the solid just dissolves.

Spectroscopy

At the option of the instructor, obtain the infrared spectrum using the dry-film method (see Technique 25, Section 25.4) or by an alternative method (see Technique 25, Section 25.5A). You should observe absorbances at 1734 and 1660 cm^{-1} for the ester carbonyl and enone groups, respectively. Compare your spectrum to that shown in Figure 1. Your instructor may also want you to determine the ^1H and ^{13}C spectra. The spectra may be run in CDCl_3 or $\text{DMSO}-d_6$. The ^1H spectrum (500 MHz CDCl_3) is shown in Figure 2. Assignments have been made on the spectrum using data from a paper by Delaude, Grandjean, and Noels (see references below.) No attempt has been made to analyze the phenyl resonances, other than to show the integral value (10 H) for the two monosubstituted benzene rings. For reference, the ^{13}C spectrum (75 MHz, CDCl_3) shows 17 peaks: 14.1, 36.3, 44.3, 59.8, 61.1, 124.3, 126.4, 127.5, 127.7, 129.0, 129.1, 130.7, 137.9, 141.2, 158.8, 169.5, and 194.3.

REFERENCES

- García-Raso, A.; García-Raso, J.; Campaner, B.; Mestres, R.; Sinisterra, J. V. An Improved Procedure for the Michael Reaction of Chalcones. *Synthesis* **1982**, 1037.
- García-Raso, A.; García-Raso, J.; Sinisterra, J. V.; Mestres, R. Michael Addition and Aldol Condensation: A Simple Teaching Model for Organic Laboratory. *J. Chem. Educ.* **1986**, 63, 443.
- Delaude, L.; Grandjean, J.; Noels, A. F. The Step-by-Step Robinson Annulation of Chalcone and Ethyl Acetoacetate. *J. Chem. Educ.* **2006**, 83, 1225–1228, and supplementary materials submitted with this article.

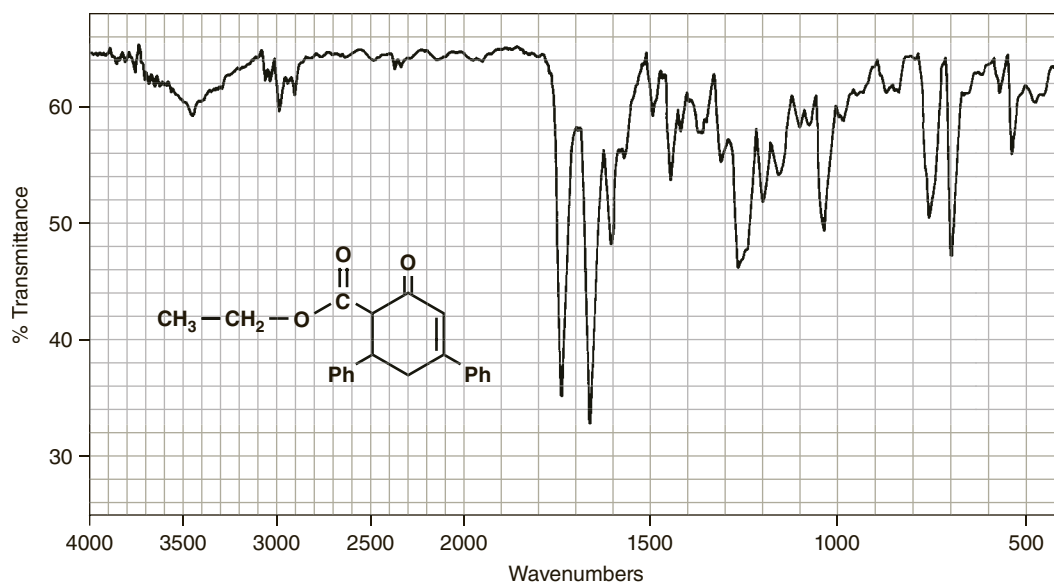


Figure 1
Infrared spectrum of 6-ethoxycarbonyl-3,5-diphenyl-2-cyclohexenone, KBr.

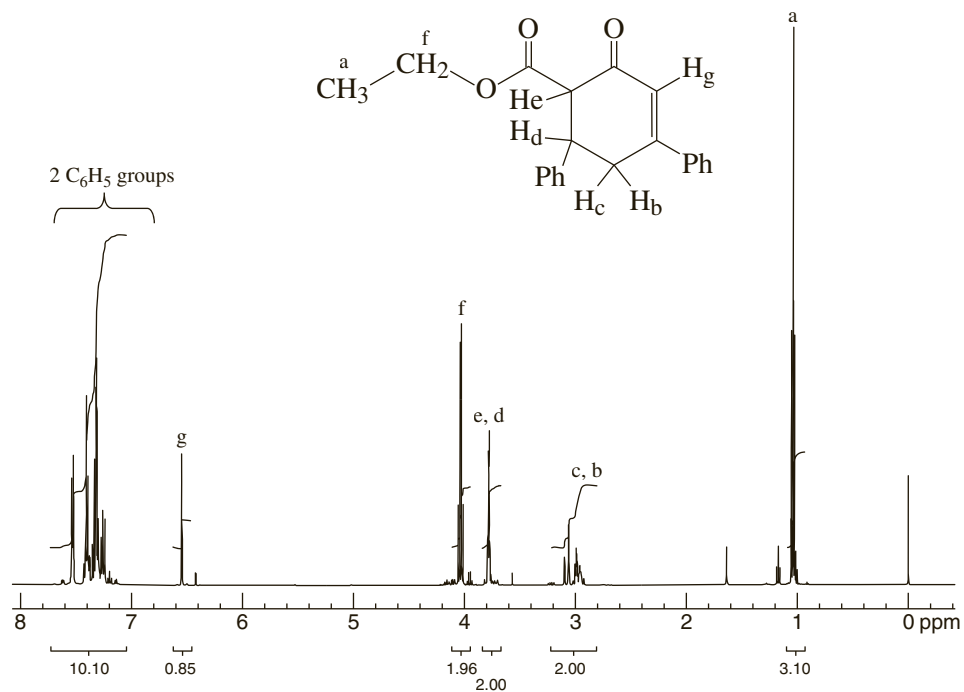


Figure 2

500 MHz ^1H NMR spectrum of 6-ethoxycarbonyl-3,5-diphenyl-2-cyclohexenone, CDCl_3 . Integral values for each of the patterns are inserted under the peaks to assign the number of protons in each pattern. Protons H_d and H_e overlap at 3.8 ppm in CDCl_3 integrating for 2H. In $\text{DMSO}-d_6$, the protons H_d and H_e are totally resolved and appear individually at 3.6 and 4.1 ppm, respectively. The other protons appear at nearly the same values in both solvents. Small impurity peaks appearing in the spectrum can be ignored.

QUESTIONS

1. Why was it possible to separate the product from sodium hydroxide using acetone?
2. The white solid that remains in the centrifuge tube after acetone extraction fizzes when hydrochloric acid is added, suggesting that sodium carbonate is present. How did this substance form? Give a balanced equation for its formation. Also give an equation for the reaction of sodium carbonate with hydrochloric acid.
3. Draw a mechanism for each of the three steps in the preparation of the 6-ethoxycarbonyl-3,5-diphenyl-2-cyclohexenone. You may assume that sodium hydroxide functions as a base and ethanol serves as a proton source.
4. Indicate how you could synthesize *trans*-chalcone. (*Hint*: Experiment 39)

Preparation of Triphenylpyridine

Green Chemistry

Aldol condensation reaction


Michael addition reaction

Crystallization

Solventless reaction

This experiment is another demonstration of a series of synthetic reactions, specifically aldol condensation followed by Michael addition, which was illustrated in Experiment 41. In this case, however, the method is designed to follow the principles of Green Chemistry. By contrast, the procedure in Experiment 41 includes the use of organic solvents (ethanol and acetone).

This experiment incorporates an aldol condensation reaction, followed by a Michael condensation, to provide a product with an interesting structure. The “green” feature in this experiment is that the entire reaction sequence is conducted without the use of any solvent at all. Avoiding the use of solvents altogether is in accord with Principle 5 of the Twelve Principles of Green Chemistry (see the essay “Green Chemistry” that precedes Experiment 29): The use of auxiliary substances (solvents, separation agents, etc.) should be avoided whenever possible and innocuous when used.

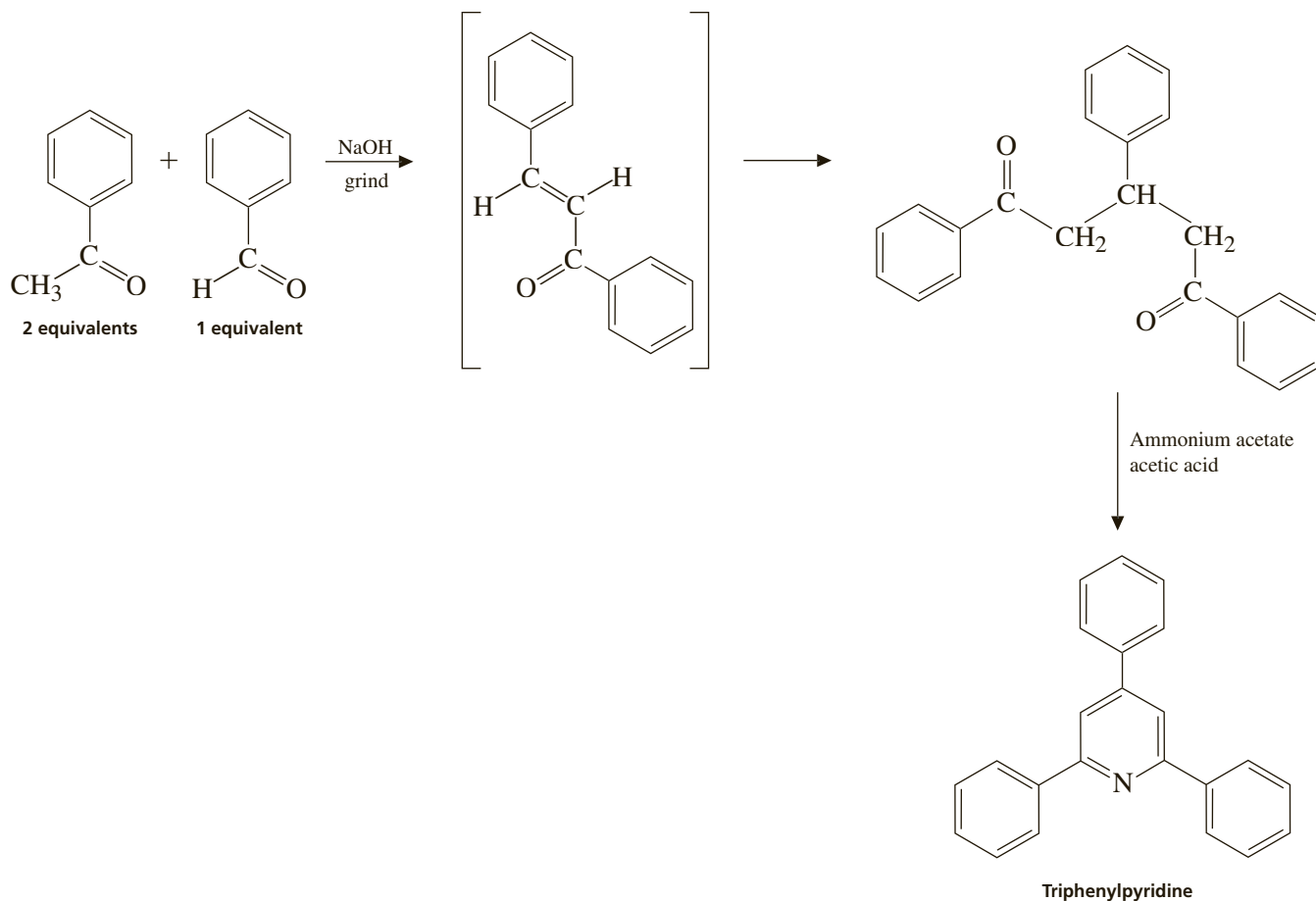

SUGGESTED WASTE DISPOSAL

All aqueous waste can be disposed of in a waste container designated for non-halogenated aqueous waste. The mortar and pestles should be rinsed with acetone, and this waste should be placed in a container intended for organic waste.

NOTES TO THE INSTRUCTOR

This procedure works best if the sodium hydroxide pellets are fresh. The quality of the product and the ease with which students will be able to grind the reagents in the mortar and pestle will be improved. It is also recommended that students work in pairs for this procedure, in order to share the workload of the lengthy period of grinding.

SAFETY PRECAUTIONS

Sodium hydroxide pellets are corrosive; they should be handled with care. Gloves should be worn during the first part of this reaction.

PROCEDURE

Part A. Michael-Aldol Condensation Reaction

To a clean dry mortar and pestle, add one sodium hydroxide pellet (0.075 to 0.095 g) and grind it to a powder. Add 0.24 g acetophenone and grind the mixture until it is homogeneous. Then add 0.11 g benzaldehyde and continue to grind.

The mixture will go through several stages, through intermediates that resemble a sticky paste, until it becomes a solid. Expect to grind (mix) the mixture for 15 minutes. Grind it thoroughly. If necessary, a metal spatula can be used to scrape the product from the sides of the mortar so the mixture can continue to be ground. Work in pairs in order to share the grinding chore to ensure that the grinding has been thorough and has continued for the full 15 minutes. Also, letting the sample stand for 15 to 20 minutes will give it time to harden. When the mixture becomes too difficult to mix, it will usually harden significantly with time, and then it can be broken up and ground. Just give it your best effort for 15 minutes, then let the mixture stand for 20 minutes. The reaction mixture must be ground very well for the full 15 minutes, but it will be messy in the beginning; mix more gently at first until it becomes more solidified or a powder, then start grinding more forcefully.

Part B. Synthesis of Triphenylpyridine

Add 0.15 g of ammonium acetate to a 25-mL round-bottom flask equipped with a stir bar. Measure 10-mL glacial acetic acid and carefully add it to the round-bottom flask. Stir this mixture for 5 minutes. Prepare a water-cooled condenser, and after transferring the product from Part A to the suspension in the round-bottom flask, connect the condenser to the flask. Heat the mixture to boiling, and allow the mixture to reflux for 2 hours. Cool the reaction mixture to room temperature with the condenser still attached. When the glassware is cool, add 10 mL of water, remove the flask from the condenser, place the flask in a labeled beaker, and place the beaker in the freezer.

Isolation of Product. After preparing a Hirsch funnel for vacuum filtration, draw 1 or 2 mL of water through the funnel to ensure a proper seal between the filter paper and funnel. Then, vacuum filter the precipitate from the round-bottom flask. Rinse the flask three times with 1-mL portions of water and also pass these portions through the vacuum filter. Transfer the product to a 25-mL Erlenmeyer flask and add 10 mL of a 5% solution of sodium bicarbonate. Swirl the mixture for 5 minutes. Transfer the product carefully, because wet filter paper tears very easily. Vacuum filter again, and then rinse the isolated precipitate twice with 1-mL portions of water.

Allow the product to stand under vacuum for 10 minutes to dry it more completely and transfer it to a watch glass to dry. Recrystallize the product from ethyl acetate. Weigh the dry triphenylpyridine and calculate the percentage yield. Determine the melting point of the product (literature value = 137°C to 138°C). Determine the proton and ¹³C NMR spectra of the product and include them and the interpretations in your laboratory report.

REFERENCES

- Palleros, D. R. Solvent-Free Synthesis of Chalcones. *J. Chem. Educ.* **2004**, *81*, 1345–1347.
- Cave, G. W. V.; Raston, C. L. Efficient Synthesis of Pyridines *via* a Sequential Solventless Aldol Condensation and Michael Addition. *J. Chem. Soc. Perkin Trans. I* **2001**, *24*, 3258–3264.

The Wittig Reaction: Preparation of 1,4-Diphenyl-1,3-butadiene

Wittig reaction

Working with sodium ethoxide

Thin-layer chromatography


UV/NMR spectroscopy

Green Chemistry


Solventless Wittig reactions

The Wittig reaction is often used to form alkenes from carbonyl compounds. In this experiment, the isomeric dienes *cis, trans*-, and *trans, trans*-1,4-diphenyl-1,3-butadiene will be formed from cinnamaldehyde and benzyltriphenyl-phosphonium chloride.




Two procedures are provided for preparing *trans, trans*-1,4-diphenyl-1,3-butadiene. In Experiment 43B, the reaction uses sodium ethoxide in ethanol solvent as the base, whereas in Experiment 43C, a green chemistry alternative is provided whereby potassium phosphate is employed as the base that is conducted without any solvent. The mechanism of the Wittig reaction in the presence of either sodium ethoxide or potassium phosphate is essentially identical. Sodium ethoxide is shown as the base in the mechanism that follows. In Experiment 43A, an optional procedure is provided for preparing one of the starting materials for the Wittig reaction.

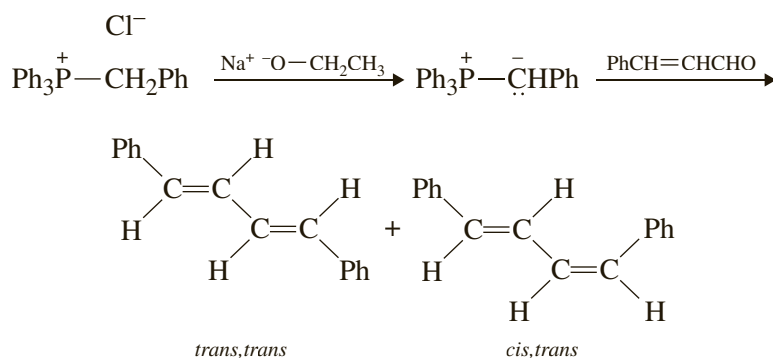




The reaction is carried out in two steps. First, the phosphonium salt is formed by the reaction of triphenylphosphine with benzyl chloride in Experiment 43A. The reaction is a simple nucleophilic displacement of chloride ion by triphenylphosphine. The salt that is formed is called the “Wittig reagent” or “Wittig salt.”

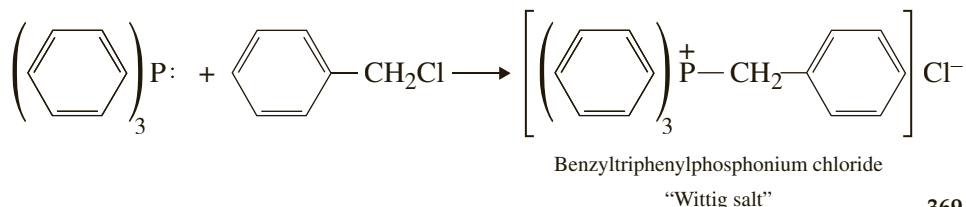

When treated with base, the Wittig salt forms an **ylide**. An ylide is a species having adjacent atoms oppositely charged. The ylide is stabilized due to the ability of phosphorus to accept more than eight electrons in its valence shell. Phosphorus uses its 3d orbitals to form the overlap with the 2p orbital of carbon that is necessary for resonance stabilization. Resonance stabilizes the carbanion.

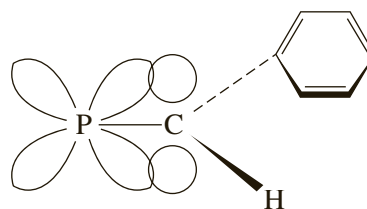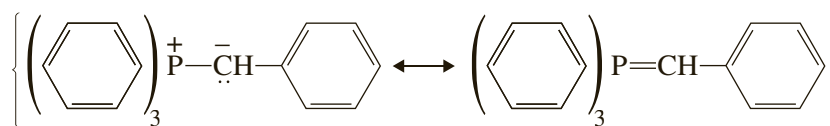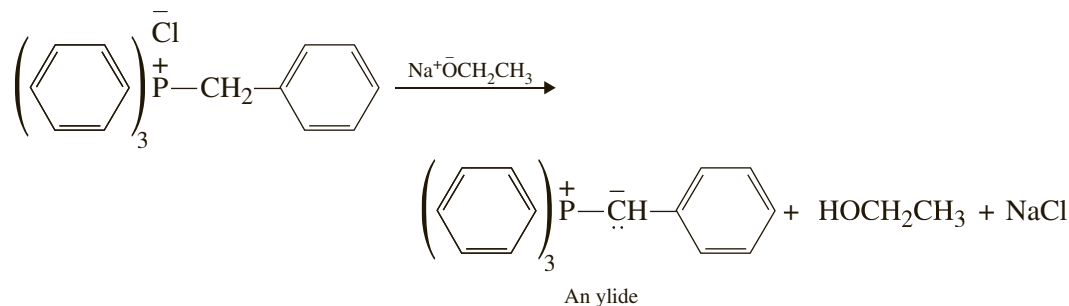

The ylide is a carbanion that acts as a nucleophile, and it adds to the carbonyl group in the first step of the mechanism. Following the initial nucleophilic addition, a remarkable sequence of events occurs, as outlined in the following mechanism:

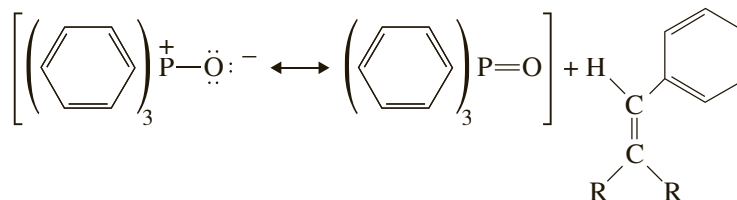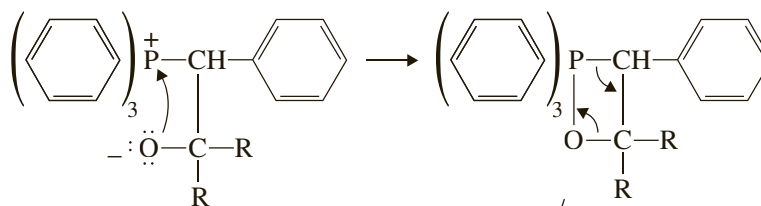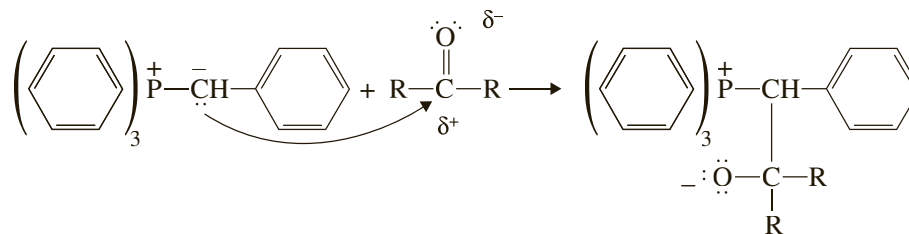


Triphenylphosphine oxide

An alkene

The addition intermediate, formed from the ylide and the carbonyl compound, cyclizes to form a four-membered-ring intermediate. This new intermediate is unstable and fragments into an alkene and triphenylphosphine oxide. Notice that the ring breaks open differently from the way it was formed. The driving force for this ring-opening process is the formation of a very stable substance, triphenylphosphine oxide. A large decrease in potential energy is achieved on the formation of this thermodynamically stable compound.

In this experiment, cinnamaldehyde is used as the carbonyl compound and yields mainly the *trans,trans*-1,4-diphenyl-1,3-butadiene, which is obtained as a solid. The *cis,trans* isomer is formed in smaller amounts, but it is a liquid that is not isolated in this experiment. The *trans,trans* isomer is the more stable isomer and is formed preferentially.

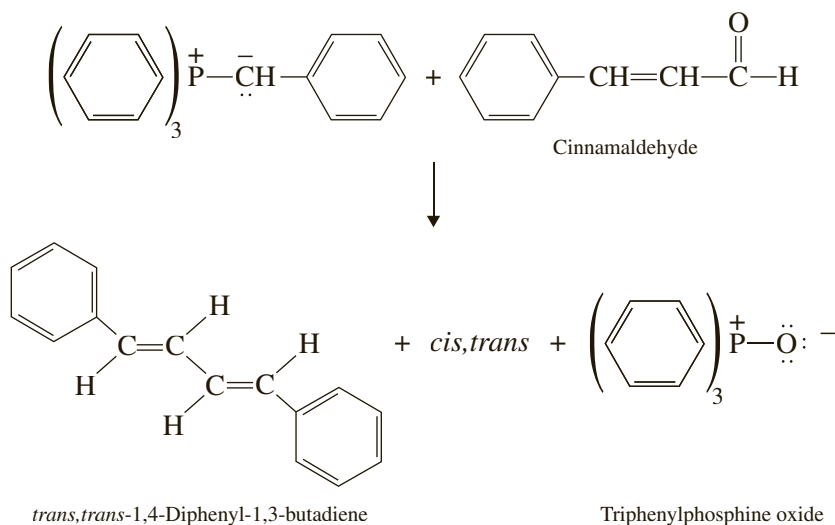

REQUIRED READING

Review: Technique 8 Section 8.3
Technique 20

SPECIAL INSTRUCTIONS

Your instructor may ask you to prepare 1,4-diphenyl-1,3-butadiene starting with commercially available benzyltriphenylphosphonium chloride. If so, start with Experiment 43B or 43C. The sodium ethoxide solution used in Experiment 43B must be kept tightly stoppered when not in use because it reacts readily with atmospheric water. Fresh cinnamaldehyde must be used in this experiment. Old cinnamaldehyde should be checked by infrared spectroscopy to be certain that it does not contain any cinnamic acid.

If your instructor asks you to prepare benzyltriphenylphosphonium chloride in Experiment 43A, you can conduct another experiment concurrently during the 1.5-hour reflux period. Triphenylphosphine is rather toxic.

Be careful not to inhale the dust. Benzyl chloride is a skin irritant and a lachrymator. It should be handled in the hood with care.

SUGGESTED WASTE DISPOSAL

If you conducted Experiment 43A, place the wastes in the nonhalogenated waste container. For Experiment 43B and 43C, dispose of organic wastes in the nonhalogenated waste container. Place the aqueous waste into the bottle designated for aqueous wastes.

43A EXPERIMENT 43A (OPTIONAL)

Benzyltriphenylphosphonium Chloride (Wittig Salt)

Place 0.550 g of triphenylphosphine ($MW = 262.3$) into a 5-mL conical vial. In a hood, transfer 0.36 mL of benzyl chloride ($MW = 126.6$, $d = 1.10$ g/mL) to the vial and add 2.0 mL of xylenes (mixture of *o*-, *m*-, and *p*-isomers).

CAUTION



Benzyl chloride is a lachrymator, a tear-producing substance.

Add a magnetic spin vane to the conical vial and attach a water-cooled condenser. Boil the mixture using an aluminum block at about 165°C for at least 1.5 hours. An increased yield may be expected when the mixture is heated longer. In fact, you may begin heating the mixture before the temperature has reached the values given, but do not include this time in the 1.5-hour reaction period. The solution will be homogeneous at first, and then the Wittig salt will begin to precipitate. Maintain the stirring during the entire heating period or bumping may occur. Remove the apparatus from the aluminum block and allow it to cool for a few minutes. Remove the vial and cool it thoroughly in an ice bath for about 5 minutes.

Collect the Wittig salt by vacuum filtration using a Hirsch funnel. Use three 1-mL portions of cold petroleum ether (bp 60–90°C) to aid the transfer and to wash the crystals free of the xylene solvent. Dry the crystals, weigh them, and calculate the percentage yield of the Wittig salt. At the option of the instructor, obtain the proton NMR spectrum of the salt in $CDCl_3$. The methylene group appears as a doublet ($J = 14$ Hz) at 5.5 ppm because of 1H - ^{31}P coupling.

43B EXPERIMENT 43B

Preparation of 1,4-Diphenyl-1,3-Butadiene Using Sodium Ethoxide to Generate the Ylide

In the following operations, cap the 5-mL conical vial whenever possible to avoid contact with moisture from the atmosphere. If you prepared your own benzyltriphenylphosphonium chloride in Experiment 43A, you may need to supplement your yield in this part of the experiment.

Preparation of the Ylide

Place 0.480 g of benzyltriphenylphosphonium chloride ($MW = 388.9$) in a dry 5-mL conical vial. Add a magnetic spin vane. Transfer 2.0 mL of absolute (anhydrous) ethanol to the vial and stir the mixture to dissolve the phosphonium salt (Wittig salt). Add 0.75 mL of sodium ethoxide solution¹ to the vial using a dry pipette while stirring continuously. Cap the vial and stir this mixture for 15 minutes. During this period, the cloudy solution acquires the characteristic yellow color of the ylide.

Reaction of the Ylide with Cinnamaldehyde

Measure 0.15 mL of pure cinnamaldehyde ($MW = 132.2$, $d = 1.11$ g/mL) and place it in another small conical vial. Add 0.50 mL of absolute ethanol to the cinnamaldehyde. Cap the vial until it is needed. After the 15-minute period, use a Pasteur pipette to mix the cinnamaldehyde with the ethanol and add this solution to the ylide in the reaction vial. A color change should be observed as the ylide reacts with the aldehyde and the product precipitates. Stir the mixture for 10 minutes.

Separation of the Isomers of 1,4-Diphenyl-1,3-Butadiene

Cool the vial thoroughly in an ice-water bath (10 min), stir the mixture with a spatula, and transfer the material from the vial to a Hirsch funnel under vacuum. Use two 1-mL portions of ice-cold absolute ethanol to aid the transfer and to rinse the product. Dry the crystalline *trans,trans*-1,4-diphenyl-1,3-butadiene by drawing air through the solid. The product has a small amount of sodium chloride that is removed as described in the next paragraph. The cloudy material in the filter flask contains triphenylphosphine oxide, the *cis,trans*-isomer, and some *trans,trans* product. Pour the filtrate into a beaker and save it for the thin-layer chromatography experiment described in the next section.

Remove the *trans,trans*-1,4-diphenyl-1,3-butadiene from the filter paper, place the solid in a 10-mL beaker, and add 3 mL of water. Stir the mixture and filter it on a Hirsch funnel, under vacuum, to collect the nearly colorless crystalline *trans,trans* product. Use about 1 mL of water to aid the transfer. Allow the solid to dry thoroughly.

Thin-Layer Chromatography

Use thin-layer chromatography to analyze the filtrate that you saved in the previous section. This mixture must be analyzed as soon as possible so that the *cis,trans* isomer will not be photochemically converted to the *trans,trans* compound. Use a 2×8 -cm silica gel TLC plate that has a fluorescent indicator (Eastman Chromatogram Sheet, No. EK 1224294). At one position on the TLC plate, spot the filtrate, as is, without dilution. Dissolve a few crystals of the *trans,trans*-1,4-diphenyl-1,

¹ This reagent is prepared in advance by the instructor. Carefully dry a 250-mL Erlenmeyer flask and insert a drying tube filled with calcium chloride into a one-hole rubber stopper. Obtain a large piece of sodium, clean it by cutting off the oxidized surface, weigh out a 2.30-g piece, cut it into 20 smaller pieces, and store it under xylene. Using tweezers, remove each piece, wipe off the xylene, and add the sodium slowly over a period of about 30 minutes to 40 mL of absolute (anhydrous) ethanol in the 250-mL Erlenmeyer flask. After the addition of each piece, replace the stopper. The ethanol will warm as the sodium reacts, but do not cool the flask. After the sodium has been added, warm the solution and shake it gently until all the sodium reacts. Cool the sodium ethoxide solution to room temperature. This reagent may be prepared in advance of the laboratory period, but it must be stored in a refrigerator between laboratory periods. When it is stored in a refrigerator, it may be kept for about 3 days. Before using this reagent, bring it to room temperature and swirl it gently in order to redissolve any precipitated sodium ethoxide. Keep the flask stoppered between each use.

3-butadiene in a few drops of acetone and spot it at another position on the plate. Use hexane as a solvent to develop (run) the plate.

Visualize the spots with a UV lamp using both the long and short wavelength settings. The order of increasing R_f values is as follows: triphenylphosphine oxide, *trans,trans*-diene, and *cis,trans*-diene. It is easy to identify the spot for the *trans,trans* isomer because it fluoresces brilliantly. What conclusion can you make about the contents of the filtrate and the purity of the *trans,trans* product? Report the results that you obtain, including R_f values and the appearance of the spots under illumination. Discard the filtrate in the container designated for nonhalogenated waste.

Yield Calculation and Melting-Point Determination

When the *trans,trans*-1,4-diphenyl-1,3-butadiene is dry, determine the melting point (literature, 151°C). Weigh the solid and determine the percentage yield. If the melting point is below 145°C, recrystallize a portion of the compound from hot 95% ethanol (20 mg/1.3 mL ethanol) in a Craig tube. Redetermine the melting point.

43C EXPERIMENT 43C

Preparation of 1,4-Diphenyl-1,3-Butadiene Using Potassium Phosphate to Generate the Ylide


Experiment 43C provides an alternative green chemistry method for preparing 1,4-diphenyl-1,3-butadiene by the Wittig reaction. No solvent is used in this experiment. Instead, the starting materials are ground together with potassium phosphate in a mortar and pestle. This experiment will demonstrate to students a more environmentally friendly method for carrying out a reaction that might be performed on a larger scale in industry.

The reaction will be accomplished by grinding cinnamaldehyde with benzyltriphenylphosphonium chloride and potassium phosphate (tribasic, K_3PO_4). This is done using a mortar and pestle. TLC will be used to analyze the crystallized *trans,trans*-1,4-diphenyl-1,3-butadiene product, as well as the filtrate from the crystallization procedure that contains both the *cis,trans* and *trans,trans*-1,4-diphenyl-1,3-butadiene isomers.

Reaction

Using an analytical balance, weigh out 309 mg of benzyltriphenylphosphonium chloride and 656 mg of potassium phosphate (tribasic, K_3PO_4) and place the solids into a clean and dry 6-cm (inside diameter) porcelain mortar with a pour lip. Using an automatic pipette, measure and add 100 μ L of cinnamaldehyde to the mixture in the mortar. Grind the mixture together for a total of 20 minutes. It is much easier to use a pestle that is long enough to grip securely in one's hand, thus saving one's fingers from getting sore or tired. At the beginning of the grinding operation, the mixture will act like putty and will have a definite yellow color. After a few minutes of grinding, the mixture starts to turn into a thick paste that adheres to the inside of the mortar and the edges of the pestle. Bend the end of a spatula as shown in the figure. This bent spatula is useful for scraping the material off of the inside of the mortar and

pestle and directing the mass into the center of the mortar. Repeat the scraping operation after every 1 to 2 minutes of grinding. Include that time in the total of 20 minutes of grinding time.

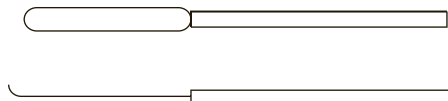

Isolation of Crude 1,4-Diphenyl-1,3-Butadiene

After 20 minutes, add a few milliliters of de-ionized water to the material in the mortar. Scrape the mortar and pestle a final time to loosen all of the product from the mortar. Pour the mixture into a Hirsch funnel inserted into a filter flask under vacuum. Use a squirt bottle with de-ionized water to transfer any remaining off-white product into the Hirsch funnel. Discard the filtrate that contains potassium phosphate and some triphenylphosphine oxide. The off-white solid consists mainly of the *trans,trans* isomer, but some of the *cis,trans* isomer will be present, as well.

Crystallization

Purify the off-white solid by crystallization from absolute ethanol in a small test tube using the standard technique of adding hot solvent until the solid dissolves. A small amount of impurity might not dissolve. If this is the case, use a Pasteur pipette to *rapidly* remove the hot solution away from the impurity and transfer the hot solution to another test tube. Cork the test tube and place it in a warm 25-mL Erlenmeyer flask. Allow the solution to cool slowly. Once the test tube has cooled and crystals have formed, place the test tube in an ice bath for at least 10 minutes to complete the crystallization process. Place 2 mL of absolute ethanol in another test tube and cool the solvent in the ice bath (this solvent will be used to aid the transfer of the product). Loosen the crystals in the test tube with a microspatula and pour the contents of the test tube into a Hirsch funnel under vacuum. Remove the remaining crystals from the test tube using the chilled ethanol and a spatula. Dry the colorless crystalline (plates) of *trans,trans*-1,4-diphenyl-1,3-butadiene on the Hirsch funnel for about 5 minutes to completely dry them. Save the filtrate from the crystallization for analysis by thin-layer chromatography. The *cis,trans*-1,4-diphenyl-1,3-butadiene, which is also formed in the Wittig reaction, is a liquid, and crystallization effectively removes the isomer from the solid *trans,trans* product.

Yield Calculation and Melting Point Determination

Weigh the purified *trans,trans* product and calculate the percentage yield. Determine the melting point of the product (literature, 151°C).

Thin-Layer Chromatography


Following the procedure in Experiment 43B, analyze the filtrate from the crystallization and the purified solid product by thin-layer chromatography. Develop the plate with hexane. This solvent will separate the *cis,trans*-diene from the *trans,trans* isomer. The order of increasing R_f values is as follows: triphenylphosphine oxide, *trans,trans*-diene, and *cis-trans*-diene. Triphenylphosphine oxide is so polar that the R_f value will be nearly zero. After developing the plate in hexane, as indicated in Experiment 43B, use the short and long wavelength settings with a UV lamp to visualize the spots. Calculate the R_f values and record them in your notebook.

Spectroscopy (Optional)

Prepare an NMR sample by dissolving at least 20 mg of crystallized product in about 1 mL of CDCl_3 in a small test tube. Transfer the solution to an NMR tube and add more solvent until the level of the solution is about 50 mm in the tube. Run the proton and carbon NMR spectra on the sample. At 300 MHz, the proton spectrum shows multiplets at 6.68 ppm and 6.95 ppm for the vinyl protons and 7.24 ppm (triplet, 2 H), 7.34 ppm (triplet, 4 H), and 7.44 ppm (doublet, 4 H) for the aromatic protons. The carbon spectrum shows peaks at 125.4, 126.6, 127.7, 128.3, 131.8, and 136.4 ppm. To determine the UV spectrum of the product, dissolve a 10-mg sample in 100 mL of hexane in a volumetric flask. Remove 10 mL of this solution and dilute it to 100 mL in another volumetric flask. This concentration should be adequate for analysis. The *trans,trans* isomer absorbs at 328 nm and possesses fine structure, whereas the *cis,trans* isomer absorbs at 313 nm and has a smooth curve.² See if your spectrum is consistent with these observations. Submit the spectral data with your laboratory report.

QUESTIONS

1. There is an additional isomer of 1,4-diphenyl-1,3-butadiene (mp 70°C), which has not been shown in this experiment. Draw the structure and name it. Why is it not produced in this experiment? (*Hint*: The cinnamaldehyde has *trans* stereochemistry.)
2. Why should the *trans,trans* isomer be the thermodynamically most stable one?
3. A lower yield of phosphonium salt is obtained in refluxing benzene than in xylene (Experiment 43A). Look up the boiling points for these solvents, and explain why the difference in boiling points might influence the yield.
4. Outline a synthesis for *cis* and *trans* stilbene (the 1,2-diphenylethenes) using the Wittig reaction.
5. The sex attractant of the female housefly (*Musca domestica*) is called **muscalure**, and its structure follows. Outline a synthesis of muscalure, using the Wittig reaction. Will your synthesis lead to the required *cis* isomer?

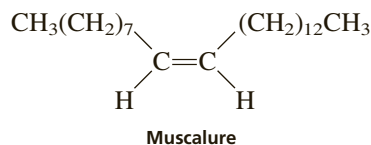

² The comparative study of the stereoisomeric 1,4-diphenyl-1,3-butadienes has been published: Pinkard, J. H., Wille, B., and Zechmeister, L. *Journal of the American Chemical Society*, 70 (1948): 1938.

Relative Reactivities of Several Aromatic Compounds


Aromatic substitution

Relative activating ability of aromatic substituents

Crystallization

When substituted benzenes undergo electrophilic aromatic substitution reactions, both the reactivity and the orientation of the electrophilic attack are affected by the nature of the original group attached to the benzene ring. Substituent groups that make the ring more reactive than benzene are called **activators**. Such groups are also said to be **ortho, para** directors because the products formed are those in which substitution occurs either ortho or para to the activating group. Various products may be formed, depending on whether substitution occurs at the ortho or para position and the number of times substitution occurs on the same molecule. Some groups may activate the benzene ring so strongly that multiple substitution consistently occurs, whereas other groups may be moderate activators, and benzene rings containing such groups may undergo only a single substitution. The purpose of this experiment is to determine the relative activating effects of several substituent groups.

In this experiment, you will study the bromination of acetanilide, aniline, and anisole:

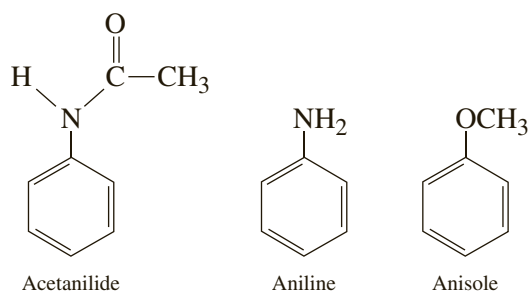


The acetamido group, $-\text{NHCOCH}_3$, the amino group, $-\text{NH}_2$, and the methoxy group, $-\text{OCH}_3$, are all activators and ortho, para directors. Each student will carry out the bromination of one of these compounds and determine its melting point. By sharing your data, you will have information on the melting points of the brominated products for acetanilide, aniline, and anisole. Using the table, it will then be possible for you to rank the three substituents in order of activating strength.


The classical method of brominating an aromatic compound is to use Br_2 and a catalyst such as FeBr_3 , which acts as a Lewis acid. The first step is the reaction between bromine and the Lewis acid:



The positive bromine ion then reacts with the benzene ring in an aromatic electrophilic substitution reaction:

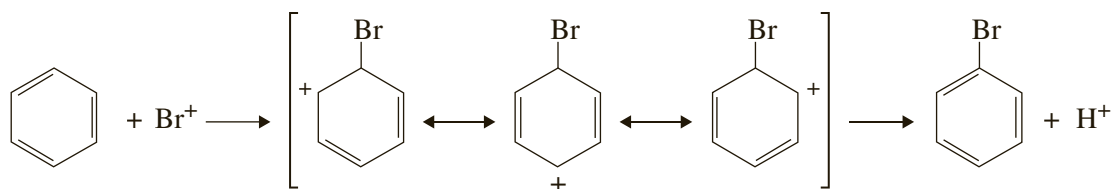

Aromatic compounds that contain activating groups can be brominated without the use of the Lewis acid catalyst because the π electrons in the benzene ring are more available and polarize the bromine molecule sufficiently to produce the required electrophile Br^+ . This is illustrated by the first step in the reaction between anisole and bromine:

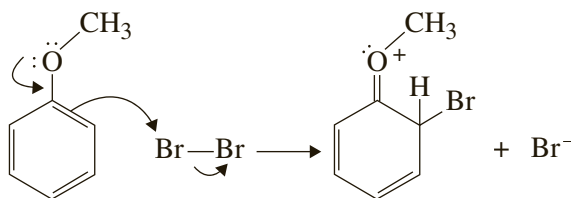

Melting points of relevant compounds

<i>Compound</i>	<i>Melting Points (°C)</i>
<i>o</i> -Bromoacetanilide	99
<i>p</i> -Bromoacetanilide	168
2,4-Dibromoacetanilide	145
2,6-Dibromoacetanilide	208
2,4,6-Tribromoacetanilide	232
<i>o</i> -Bromoaniline	32
<i>p</i> -Bromoaniline	66
2,4-Dibromoaniline	80
2,6-Dibromoaniline	87
2,4,6-Tribromoaniline	122
<i>o</i> -Bromoanisole	3
<i>p</i> -Bromoanisole	13
2,4-Dibromoanisole	60
2,6-Dibromoanisole	13
2,4,6-Tribromoanisole	87

In this experiment, the brominating mixture consists of bromine, hydrobromic acid HBr, and acetic acid. The presence of bromide ion from the hydrobromic acid helps to solubilize the bromine and increase the concentration of the electrophile.

REQUIRED READING

Review: Technique 11

You should review the chapters in your lecture textbook that deal with electrophilic aromatic substitution. Pay special attention to halogenation reactions and the effect of activating groups.

SPECIAL INSTRUCTIONS

Bromine is a skin irritant, and its vapors cause severe irritation to the respiratory tract. It will also oxidize many pieces of jewelry. Hydrobromic acid may cause skin or eye irritation. Aniline is highly toxic and a suspected teratogen. All bromoanilines are toxic. This experiment should be carried out in a fume hood or in a well-ventilated laboratory.

Each person will carry out the bromination of only one of the aromatic compounds according to your instructor's directions. The procedures are identical except for the initial compound used and the final recrystallization step.

NOTE TO THE INSTRUCTOR: Prepare the brominating mixture in advance.

SUGGESTED WASTE DISPOSAL

Dispose of the filtrate from the Hirsch funnel filtration of the crude product into a container specifically designated for this mixture. Place all other filtrates into the container for halogenated organic solvents.

PROCEDURE

Running the Reaction

To a tared 5-mL conical vial with a cap, add the given amount of one of the following compounds: 0.090 g of acetanilide, 0.060 mL of aniline, or 0.070 mL of anisole. Reweigh the conical vial to determine the actual weight of the aromatic compound. Add 0.5 mL of glacial acetic acid and a spin vane to the conical vial. Attach an air condenser and place the conical vial in a water bath at 23–27°C, as shown in Technique 6, Figure 6.6. Stir the mixture until the aromatic compound is completely dissolved. While the compound is dissolving, pack a drying tube loosely with glass wool. Add about 0.5 mL of 1M sodium bisulfite dropwise to the glass wool until it is moistened but not soaked. This apparatus will capture any bromine given off during the following reaction.

Under the hood, obtain 1.0 mL of the bromine/hydrobromic acid mixture¹ in a 3-mL conical vial. Place the cap on the vial before returning to your lab bench. While stirring, add all the bromine/hydrobromic acid mixture through the top of the air condenser, using a Pasteur pipette.

CAUTION



Be careful not to spill any of this mixture.

Attach the drying tube prepared above. Continue stirring the reaction mixture for 20 minutes.

¹ *Note to the Instructor:* The brominating mixture is prepared by adding 2.6 mL of bromine to 17.4 mL of 48% hydrobromic acid. This will provide enough solution for 20 students, assuming no waste of any type. This solution should be stored in the hood.

Crystallization and Isolation of Product

When the reaction is complete, transfer the mixture to a 10-mL Erlenmeyer flask containing 5 mL of water and 0.5 mL of saturated sodium bisulfite solution. Stir this mixture with a glass stirring rod until the red color of bromine disappears.² If an oil has formed, it may be necessary to stir the mixture for several minutes. Place the Erlenmeyer flask in an ice bath for 10 minutes. If the product does not solidify, scratch the bottom of the flask with a glass stirring rod to induce crystallization. It may take 10–15 minutes to induce crystallization of the brominated anisole product.³ Filter the product on a Hirsch funnel with suction and rinse with several 1-mL portions of cold water. Air-dry the product on the funnel for about 5 minutes with the vacuum on.

Recrystallization and Melting Point of Product

If you started with **aniline**, transfer the solid to a 10-mL Erlenmeyer flask and recrystallize the product from 95% ethanol (see Technique 11, Section 11.3, and Figure 11.4). Filter the crystals on a Hirsch funnel and dry them for several minutes with suction. The brominated products from either **acetanilide** or **anisole** should be crystallized using a Craig tube (Technique 11, Section 11.4 and Figure 11.6). Use 95% ethanol to crystallize the acetanilide product and hexane to crystallize the brominated anisole compound. Allow the crystals to air-dry and determine the weight and melting point.

According to the melting point and the preceding table, you should be able to identify your product. Calculate the percentage yield and submit your product, along with your report, to your instructor.

REPORT

By collecting data from other students, you should be able to determine which product was obtained from the bromination of each of the three aromatic compounds. Using this information, arrange the three substituent groups (acetamido, amino, and methoxy) in order of decreasing ability to activate the benzene ring.

REFERENCE

Zaczek, N. M.; Tyszkiewicz, R. B. Relative Activating Ability of Various Ortho, Para-Directors. *Journal of Chemical Education*, **1986**, 63, 510.

QUESTIONS

- Using resonance structures, show why the amino group is activating. Consider an attack by the electrophile E^+ at the *para* position.
- For the substituent in this experiment that was found to be least activating, explain why bromination took place at the position on the ring indicated by the experimental results.
- What other experimental techniques (including spectroscopy) might be used to identify the products in this experiment?

² If the color of bromine is still present, add a few more drops of saturated sodium bisulfite and stir the mixture for a few more minutes. The entire mixture, including liquid and solid (or oil), should be colorless.


³ If crystals fail to form after 15 minutes, it may be necessary to seed the mixture with a small crystal of product.

Nitration of Methyl Benzoate

Aromatic substitution


Crystallization

The nitration of methyl benzoate to prepare methyl *m*-nitrobenzoate is an example of an electrophilic aromatic substitution reaction, in which a proton of the aromatic ring is replaced by a nitro group:

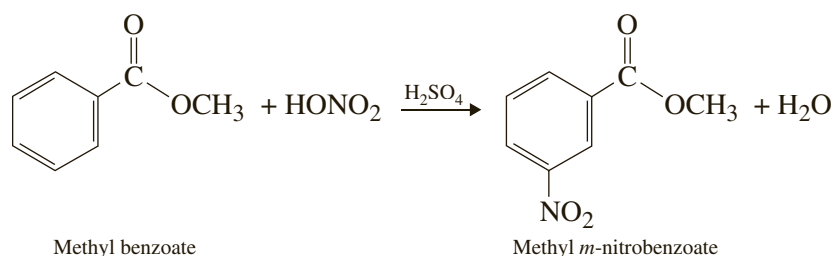

Many such aromatic substitution reactions are known to occur when an aromatic substrate is allowed to react with a suitable electrophilic reagent, and many other groups besides nitro may be introduced into the ring.

You may recall that alkenes (which are electron-rich due to an excess of electrons in the π system) can react with an electrophilic reagent. The intermediate formed is electron-deficient. It reacts with the nucleophile to complete the reaction. The overall sequence is called **electrophilic addition**. Addition of HX to cyclohexene is an example.

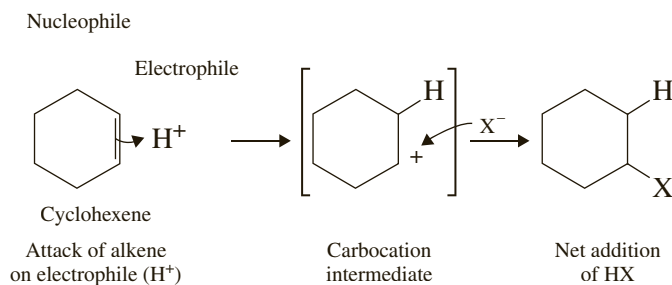

Aromatic compounds are not fundamentally different from cyclohexene. They can also react with electrophiles. However, because of resonance in the ring, the electrons of the π system are generally less available for addition reactions because an addition would mean the loss of the stabilization that resonance provides. In practice, this means that aromatic compounds react only with **powerfully electrophilic reagents**, usually at somewhat elevated temperatures.

solution. Water interferes because it is a stronger base than H_2SO_4 or HNO_3 . Temperature is also a factor in determining the extent of nitration. The higher the temperature, the greater will be the amounts of dinitration products formed in the reaction.

REQUIRED READING

Review: Technique 11 Crystallization: Purification of Solids
Technique 25 Infrared Spectroscopy, Sections 25.4 and 25.5

SPECIAL INSTRUCTIONS

It is important that the temperature of the reaction mixture be maintained below 15°C . Nitric acid and sulfuric acid, especially when mixed, are corrosive substances. Be careful not to get these acids on your skin. If you do get some of these acids on your skin, flush the affected area liberally with water.

SUGGESTED WASTE DISPOSAL

The filtrate from the Hirsch funnel filtration should be placed in the designated container.

PROCEDURE

Add 0.210 mL of methyl benzoate to a tared 3-mL conical vial and determine the actual weight of methyl benzoate. Add 0.45 mL of concentrated sulfuric acid to the methyl benzoate, along with a magnetic spin vane. Attach an air condenser to the conical vial. The purpose of the air condenser is to make it easier to hold the conical vial in place. Prepare an ice bath in a 250-mL beaker using both ice and water. Clamp the air condenser so that the conical vial is immersed in the ice bath as shown in Technique 6, Figure 6.6. (Note that in Figure 6.6 a water bath is shown rather than an ice bath.) While stirring, *very slowly* add a cool mixture of 0.15 mL of concentrated sulfuric acid and 0.15 mL of concentrated nitric acid throughout a period of about 15 minutes. The acid mixture should be added with a 9-inch Pasteur pipette through the top of the air condenser. If the addition is too fast, the formation of by-product increases rapidly, reducing the yield of the desired product.

After you have added all the acid, warm the mixture to room temperature by replacing the ice water in the 250-mL beaker with water at room temperature. Let the reaction mixture stand for 15 more minutes without stirring. Then, using a Pasteur pipette, transfer the reaction mixture to a 20-mL beaker containing 2.0 g of crushed ice. After the ice has melted, isolate the product by vacuum filtration using a Hirsch funnel and wash it with two 1.0-mL portions of cold water and then with two 0.3-mL portions of ice-cold methanol. Weigh the crude, dry product and recrystallize it from methanol using a Craig tube (see Technique 11, Section 11.4).

Determine the melting point of the product. The melting point of the recrystallized product should be 78°C . Obtain the infrared spectrum using the dry film method (Technique 25, Section 25.4) or by an alternative method. (Technique 25,


Section 25.5). Compare your infrared spectrum with the one reproduced here. Calculate the percentage yield and submit the product to the instructor in a labeled vial.

MOLECULAR MODELING (OPTIONAL)

If you are working alone, complete Part One. If you have a partner, one of you should complete Part One and the other complete Part Two. If you work with a partner, you should combine results at the end of the experiment.

Part One: Nitration of Methyl Benzoate

In this exercise, we will try to explain the observed outcome of the nitration of methyl benzoate. The major product of this reaction is methyl *m*-nitrobenzoate, where the nitro group has been added to the *meta* position of the ring. The rate-determining step of this reaction is the attack of the nitronium ion on the benzene ring. Three benzenium ion intermediates (*ortho*, *meta*, and *para*) are possible:

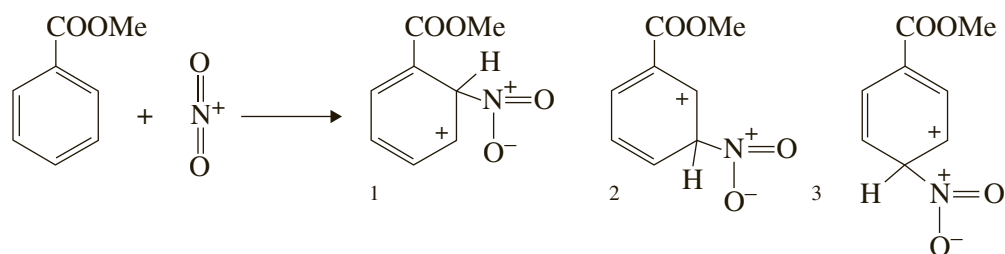

We will calculate the heats of formation for these intermediates to determine which of the three has the lowest energy. Assume that the activation energies are similar to the energies of the intermediates themselves. This is an application of the Hammond Postulate, which states that the activation energy leading to an intermediate of higher energy will be higher than the activation energy leading to an intermediate of lower energy, and vice versa. Although there are prominent exceptions, this postulate is generally true.

Make models of each of the three benzenium ion intermediates (separately), and calculate their heats of formation using an AM1-level calculation with geometry optimization. Don't forget to specify a positive charge when you submit the calculation. What do you conclude?

Now take a piece of paper and draw the resonance structures that are possible for each intermediate. Do not worry about structures involving the nitro group; only consider where the charge in the ring may be delocalized. Also note the polarity of the carbonyl group by placing a δ^+ symbol on the carbon and a δ^- symbol on the oxygen. What do you conclude from your resonance analysis?

Part Two: Nitration of Anisole


For this computation, you will analyze the three benzenium ions formed from anisole (methoxybenzene) and the nitronium ion (see Part One). Calculate the heats of formation using AM1-level calculations with geometry optimization. Don't forget to specify a positive charge. What do you conclude for anisole? How do the results compare to those for methyl benzoate?

Now take a piece of paper and draw the resonance structures that are possible for each intermediate. Do not worry about structures involving the nitro group; only consider where the charge in the ring may be delocalized. Do not forget that

the electrons on the oxygen can participate in the resonance. What do you conclude from your resonance analysis?

QUESTIONS

1. Why is methyl *m*-nitrobenzoate formed in this reaction instead of the ortho or para isomers?
2. Why does the amount of the dinitration increase at high temperatures?
3. Why is it important to add the nitric acid–sulfuric acid mixture slowly during a 15-minute period?
4. Interpret the infrared spectrum of methyl *m*-nitrobenzoate.
5. Indicate the product formed on nitration of each of the following compounds: benzene, toluene, chlorobenzene, and benzoic acid.

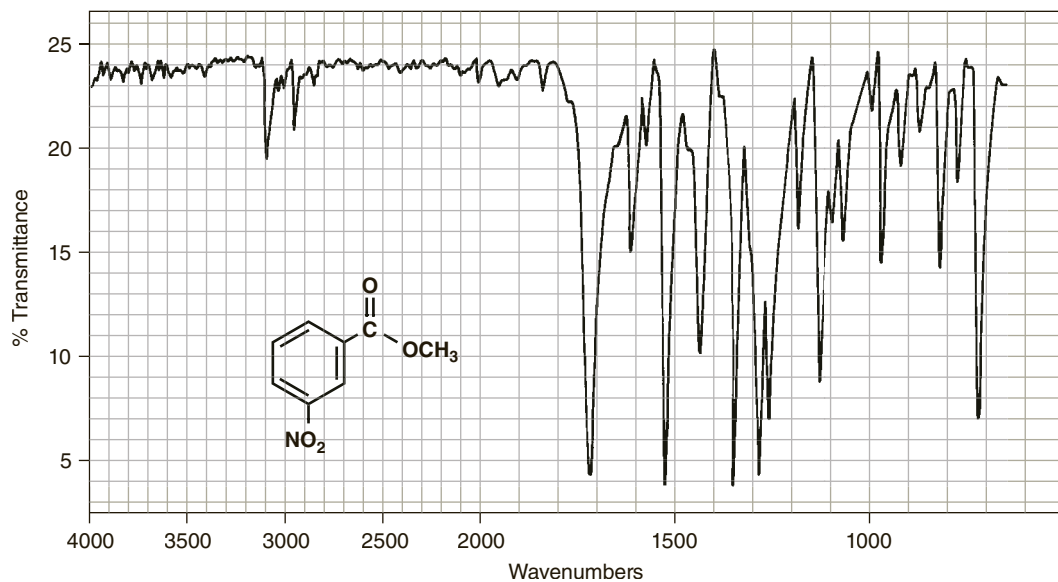

Infrared spectrum of methyl *m*-nitrobenzoate, KBr.


Synthetic Dyes

The practice of using dyes is an ancient art. Substantial evidence exists that plant dyes were known long before humans began to keep written history. Before the twentieth century, practically all dyes were obtained from natural plant or animal sources. Dyeing was a complicated and secret art, passed from one generation to the next. Dyes were extracted from plants mainly by macerating the roots, leaves, or berries in water. The extract was often boiled and then strained before use. In some cases, it was necessary to make the extraction mixture acidic or basic before the dye could be liberated from the plant tissues. Applying the dyes to cloth was also a complicated process. **Mordants** were used to fix the dye to the cloth or even to modify its color.


Madder is one of the oldest known dyes. Alexander the Great was reputed to have used the dye to trick the Persians into overconfidence during a critical battle. Using madder, a root bearing a brilliant red dye, he simulated bloodstains on the tunics of his soldiers. The Persians, seeing the apparently incapacitated Greek Army, became overconfident and much to their surprise, were overwhelmingly defeated. Through modern chemical analysis, we now know the structure of the dye found in madder root. It is called **alizarin** (see the structures below) and is very similar in structure to another ancient dye, **henna**, which has been responsible for a long line of synthetic redheads. Madder is obtained from the plant *Rubia tinctorum*. Henna is a dye prepared from the leaves of the Indian henna plant (*Lawsonia alba*) and an extract of *Acacia catechu*.

Indigo is another plant dye with a long history. This dye, obtained from the plant *Indigofera tinctoria*, has been known in Asia for more than 4000 years. By the ancient process for producing indigo, the leaves of the indigo plant are cut and illustration allowed to ferment in water. During the fermentation, **indican** (see scheme below) is extracted into the solution, and the attached glucose molecule is split off to produce **indoxyl**. The fermented mixture is transferred to large open vats, in which the liquid is beaten with bamboo sticks. During this process, the indoxyl is air-oxidized to indigo. Indigo, a strong blue dye, is insoluble in water, and it precipitates. Today, indigo is made synthetically, and its principal use is in dyeing denim to produce “blue jean” material.


Many plants yield dyes that will dye wool or silk, but there are few of these that dye cotton well. Most do not dye synthetic fibers like polyester or rayon. In addition, the natural dyes, with a few exceptions, do not cover a wide range of colors, nor do they yield “brilliant” colors. Even though some people prefer the softness of the “homespun” colors from natural dyes, the **synthetic dyes**, which give rise to deep, brilliant colors, are much in demand today. Also, synthetic dyes that will dye popular synthetic fibers can now be manufactured. Thus today, we have available an almost infinite variety of colors as well as dyes to dye any type of fabric.

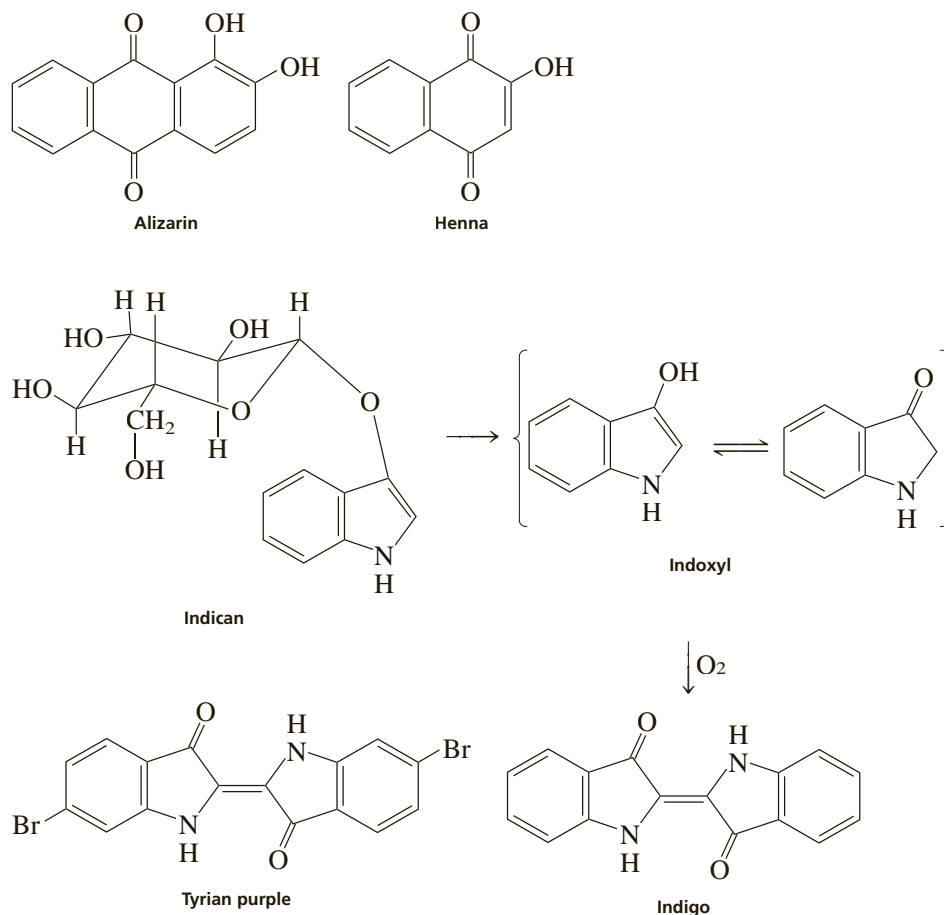

Before 1856, all dyes in use came from natural sources. However, an accidental discovery by W. H. Perkin, an English chemist, started the development of a huge synthetic dye industry, mostly in England and Germany. Perkin, then only aged 18, was trying to synthesize quinine. Structural organic chemistry was not very well developed at that time, and the chief guide to the structure of a compound was its molecular formula. Perkin thought, judging from the formulas, that it might be possible to synthesize quinine by the oxidation of allyltoluidine:

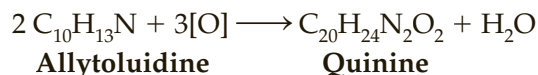


He made allyltoluidine and oxidized it with potassium dichromate. The reaction was unsuccessful, because allyltoluidine bore no structural relation to quinine. He obtained no quinine, but he did recover a reddish-brown precipitate with properties that interested him. He decided to try the reaction with a simpler base, aniline. On treating aniline sulfate with potassium dichromate, he obtained a black precipitate, which could be extracted with ethanol to give a beautiful purple solution. This purple solution subsequently proved to be a good dye for fabrics. After receiving favorable comments from dyers, Perkin resigned his post at the Royal College and went on to found the British coal tar dye industry. He became a very successful industrialist and retired at age 36 (!) to devote full time to research. The dye he synthesized became known as **mauve**. The structure of mauve was not proven until much later. From the structure (see figure) it is clear that the aniline Perkin used was not pure and that it contained the *o*-, *m*-, and *p*-toluidines also.

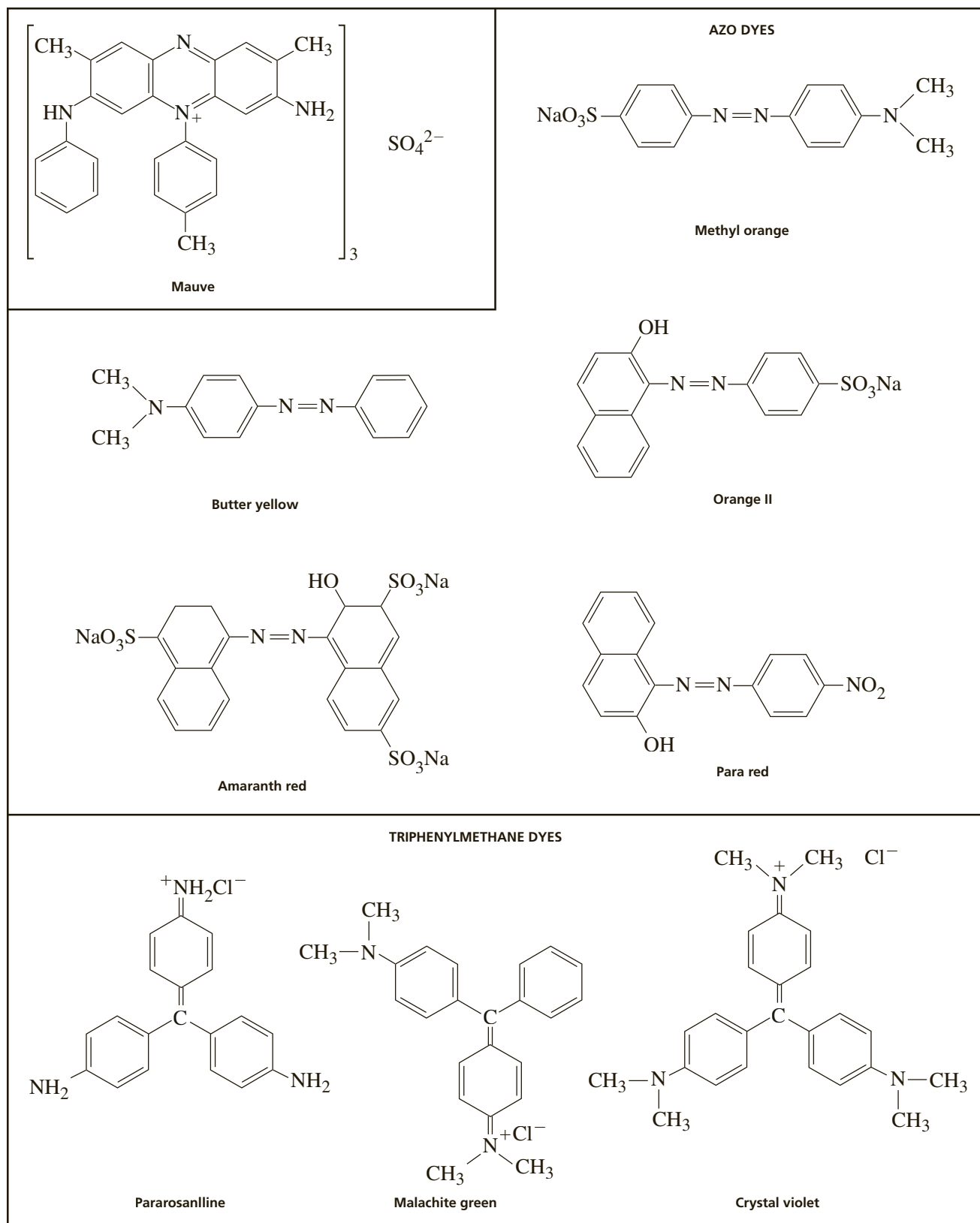


Mauve was the first synthetic, but soon (1859) the triphenylmethyl dyes pararosaniline, malachite green, and crystal violet (see figures) were discovered in France. These dyes were produced by treating mixtures of aniline, the toluidines, or both, with nitrobenzene, an oxidizing agent, and in a second step, with concentrated hydrochloric acid. The triphenylmethyl dyes were soon joined by **synthetic** alizarin (Lieberman, 1868), **synthetic** indigo (Baeyer, 1870), and the azo dyes (Griess, 1862). The azo dyes, also manufactured from aromatic amines, revolutionized the dye industry.

Baeyer's first method for producing indigo used isatin as the starting material. However, isatin was very expensive, and this method was not suitable for making large quantities of indigo. Because the commercial value of indigo was so great, much work was done to develop a less expensive method. In 1892, Baeyer and Drewson developed a method using *o*-nitrobenzaldehyde as the starting material. At the time, this method was too expensive because the synthesis of *o*-nitrobenzaldehyde was too difficult. However, *o*-nitrobenzaldehyde is now readily available, and this method will be used to synthesize indigo in Experiment 45.

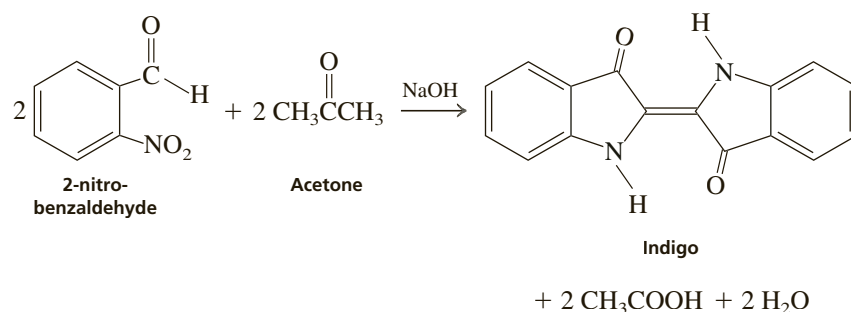

Dyes may be classified according to the chemical composition of the dyes or the way in which the dye is applied to the fabric. Indigo belongs to a category of dyes referred to as **vat dyes**. Vat dyes have a low solubility in water, and therefore, they cannot be applied directly to the fabric. When indigo undergoes reduction in base, the reduced form, leucoindigo, is soluble in water, and it is colorless. When the appropriate fabric is placed into this solution, the molecules of leucoindigo are absorbed by the pores of the fibers of the fabric and trapped in the fiber. The fabric is then removed from the "vat," rinsed thoroughly with water, and allowed to dry in the air. While drying in the presence of heat or light, the reduced form of indigo undergoes oxidation back to its insoluble form with the characteristic color of indigo.

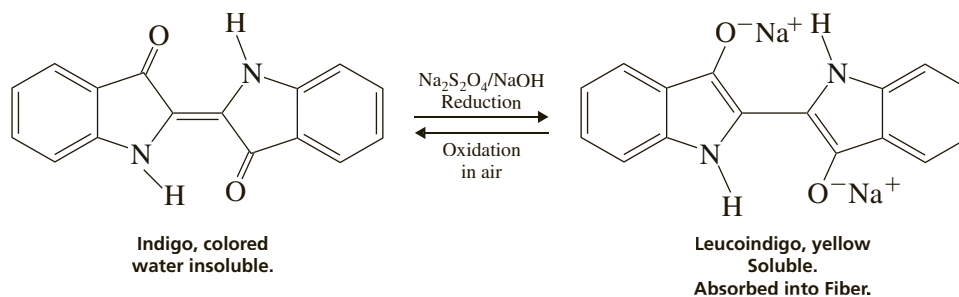


The **azo** dyes (classified according to composition) are one of the most common types of dye used today. They are used as dyes for clothing, as food dyes, and as

pigments in paints. In addition, they are used in printing inks and in certain color printing processes. Azo dyes have the basic structure




Several of these are illustrated in the figure. The unit containing the nitrogen–nitrogen bond is called an **azo group**, a strong chromophore that imparts a brilliant color to these compounds.


Producing an azo dye involves treating an aromatic amine with nitrous acid to give a **diazonium ion** intermediate. This process is called **diazotization**:

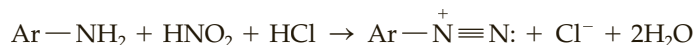

The diazonium ion is an electron-deficient (electrophilic) intermediate. A nucleophilic aromatic compound will react with the diazonium ion. The most common nucleophilic compounds are aromatic amines, and phenols. All of these types of compounds are usually more nucleophilic at a ring carbon than at either nitrogen or oxygen. This is due to resonance of the following types:

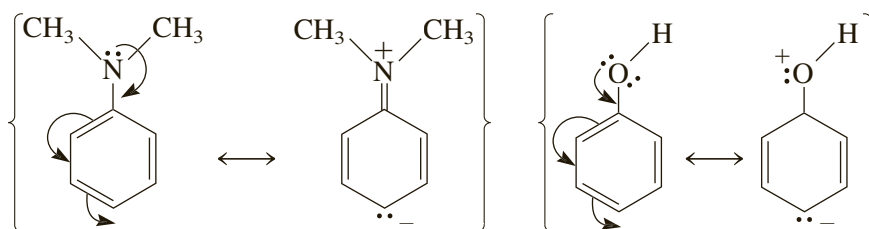

The addition of the amine or the phenol to the diazonium ion is called the **diazonium coupling** reaction, and it takes place as shown:

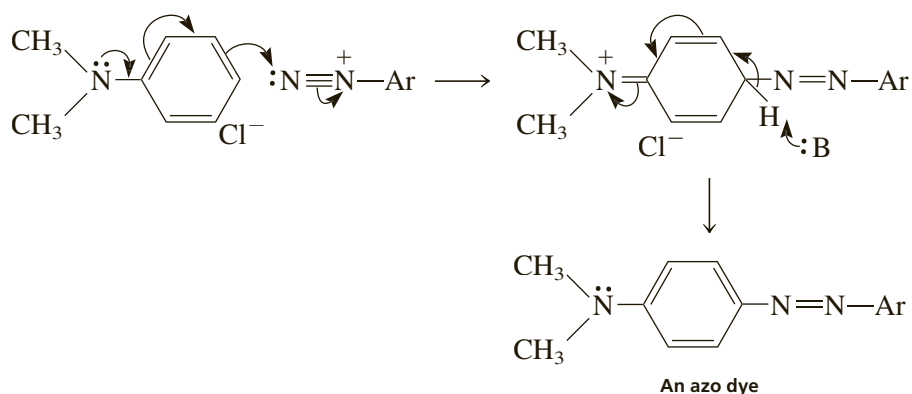

Azo dyes are both the largest and the most important group of synthetic dyes. In the formation of the azo linkage, many combinations of ArNH_2 and $\text{Ar}'\text{NH}_2$ [or $\text{Ar}'(\text{OH})$] are possible. These combinations give rise to dyes with a broad range of colors, encompassing yellows, oranges, reds, browns, and blues. The preparation of an azo dye is given in Experiment 46.

Although azo dyes are most commonly used as food colors or in dyeing fabrics, they can also be used as pigments in paints. In Experiment 48, a procedure is presented for formulating a paint using various dyes as the pigment, including methyl orange (from Experiment 46) and indigo (from Experiment 47).

The azo dyes, the triphenylmethyl dyes, and mauve are all synthesized from the anilines (aniline, *o*-, *m*-, and *p*-toluidine) and aromatic compounds (benzene, naphthalene, anthracene). All of these substances can be found in **coal tar**, a crude

material that is obtained by distilling coal. Perkin's discovery led to a multimillion-dollar industry based on coal tar, a material that was once widely regarded as a foul-smelling nuisance. Today, these same materials can be recovered from crude oil or from petroleum as by-products in the refining of gasoline. Although we no longer use coal tar, many of the dyes are still widely used.

Azo dyes account for about 60%–70% of all dyes used in food and fabrics. The toxicity of azo dyes is very low (low LD₅₀ values), and only a very small amount is required to impart the desired color. Since they are water soluble, they do not accumulate in the body. However, a few of the azo dyes break down in the body or environment to form carcinogenic products. Because of these concerns, in 2003 the European Union banned the use of 22 azo dyes containing aromatic amines. In 2010, the EPA issued an action plan to address similar concerns. This action plan requires that new products containing benzidine-based azo dyes must be submitted to the EPA for review prior to going on the market.

REFERENCES

- Abraham, E. N. *Dyes and Their Intermediates*; Pergamon Press: London, 1968.
- Allen, R. L. M. *Colour Chemistry*; Appleton-Century-Crofts: New York, 1971.
- Decelles, C. The Story of Dyes and Dyeing. *J. Chem. Educ.* **1949**, *26*, 583.
- Jaffee, H. E.; Marshall, W. J. The Origin of Color in Organic Compounds. *Chemistry* **1964**, *37* (Dec), 6.
- Juster, N. J. Color and Chemical Constitution. *J. of Chem. Edu.* **1957**, *34*, 596.
- McCoy, M. Into the Blue. *Chem. and Eng. News* **2012**, *90* (April 2), 18.
- Robinson, R. Sir William Henry Perkin: Pioneer of Chemical Industry. *J. Chem. Educ.* **1957**, *34*, 54.
- Sementsov, A. The Medical heritage from Dyes. *Chemistry* **1966**, *39* (Nov), 20.
- Shaar, B. E. Chance Favors the Prepared Mind. Part 1: Aniline Dyes. *Chemistry* **1966**, *39* (Jan), 12.
- www.food-info.net/uk/colour/azo.htm
- www.epa.gov/oppt/existingchemicals/pubs/actionplans/benzidine.html


Preparation of Methyl Orange

Diazotization

Diazonium coupling


Azo dyes

In this experiment, the azo dye **methyl orange** is prepared by the diazo coupling reaction. It is prepared from sulfanilic acid and *N,N*-dimethylaniline. The first product obtained from the coupling is the bright red acid form of methyl orange, called **helianthin**. In base, helianthin is converted to the orange sodium salt, called methyl orange. The basic scheme is shown below.

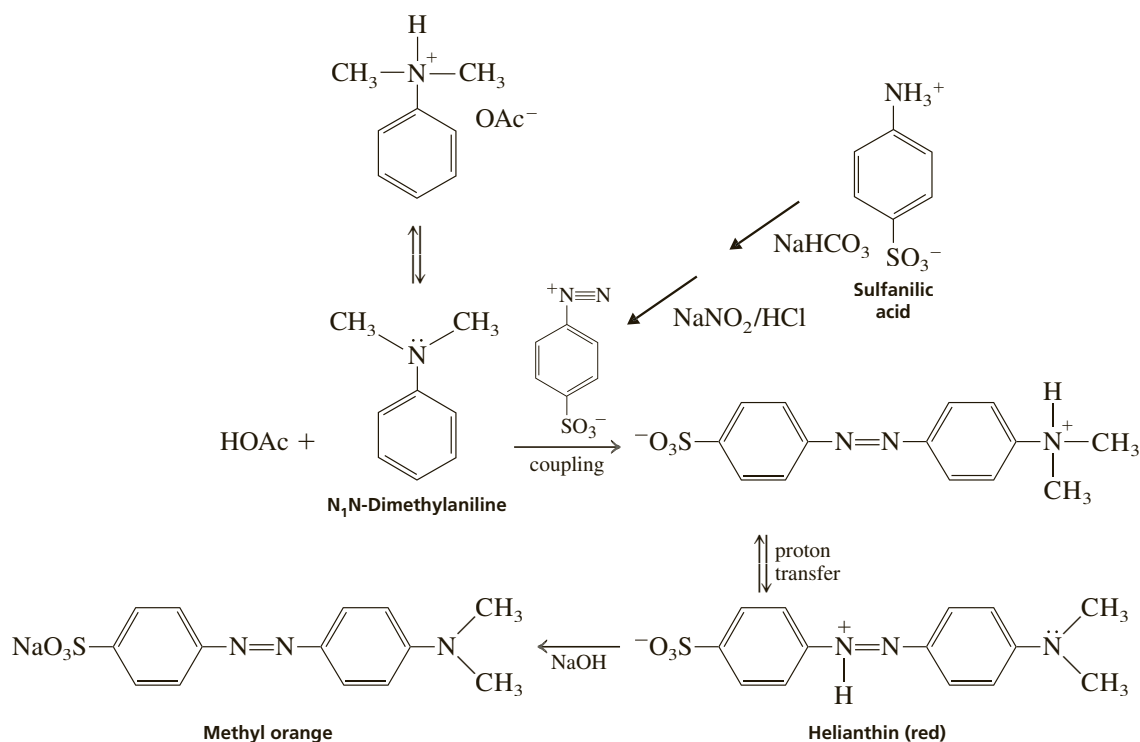

Although the sulfanilic acid is insoluble in acid solutions, it is nevertheless necessary to carry out the diazotization reaction in an acid solution. This problem can be avoided by first converting the sulfanilic acid to its soluble sodium salt by treatment with a weak base (NaHCO_3). When the basic solution is acidified, the sulfanilic acid precipitates as a fine suspension that reacts instantly with nitrous acid.


Step One. Dissolve sulfanilic acid in basic solution.

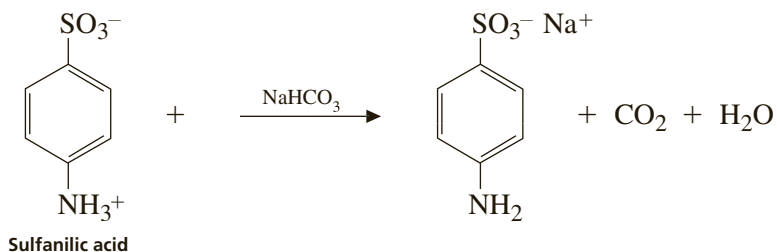

Step Two. When sodium nitrate is added and the solution is acidified, nitrous acid forms.




In the acidified solution, the sulfanilic acid is precipitated out of solution as a finely divided solid, which is immediately diazotized.

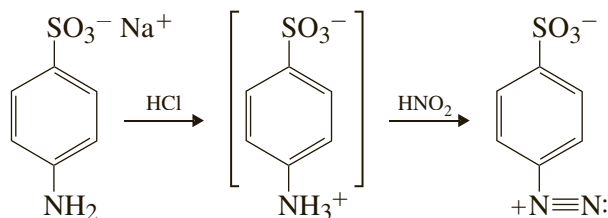


Step Three. The finely divided diazonium salt reacts immediately with added dimethylaniline, as shown in the basic scheme at the top of this experiment. The mechanism of the formation of the azo compound is shown in the essay "Synthetic Dyes."

Methyl orange is often used as an acid–base indicator. In solutions that are more basic than pH 4.4, methyl orange exists almost entirely as the **yellow** negative ion. In solutions that are more acidic than pH 3.2, it is protonated to form a **red** dipolar ion.

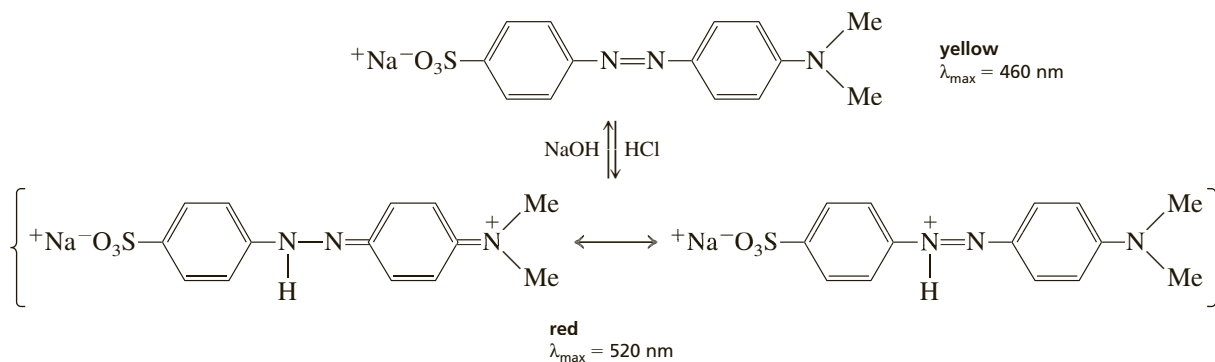

Thus, methyl orange can be used as an indicator for titrations that have their end point in the pH 3.2 to 4.4 region. The indicator is usually prepared as a 0.01% solution in water. In higher concentrations in basic solution, of course, methyl orange appears **orange**.

Azo compounds are easily reduced at the nitrogen–nitrogen double bond by reducing agents. Sodium hydrosulfite, $\text{Na}_2\text{S}_2\text{O}_4$, is often used to bleach azo compounds:

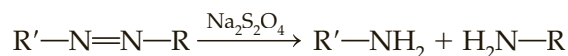

Other good reducing agents, such as stannous chloride in concentrated hydrochloric acid, will also work.

REQUIRED READING

<i>Review:</i>	Technique 8	Filtration, Section 8.3
<i>New:</i>	Essay	Synthetic Dyes

SPECIAL INSTRUCTIONS

Temperature is important in this reaction. You must keep the temperature of the reaction mixture below 10°C when the diazonium salt is formed. If the diazonium solution is stored for even a few minutes, it should be kept in an ice bath. If the temperature of the aqueous diazonium salt solution is allowed to rise, the diazonium salt will be hydrolyzed to a phenol, reducing the yield of the desired product.

Do not allow *N,N*-dimethylaniline to come in contact with your skin. When using this substance always wear rubber gloves and open the bottle in the hood. *N,N*-Dimethylaniline is a known carcinogen.

NOTES TO THE INSTRUCTOR

The methyl orange synthesized in this experiment can be used in Experiment 48, Formulation of a Paint and Art Project.

WASTE DISPOSAL

Dispense of any organic materials in the waste container for nonhalogenated organic waste. Dispose of any aqueous solutions in the waste container designated for the disposal of aqueous solutions.

PROCEDURE

Diazotized Sulfanilic Acid. Add 0.12 g of sodium bicarbonate and 10 mL of water to a 50-mL Erlenmeyer flask. Stir with a magnetic stirring bar until all of the solid is dissolved. Add 0.36 g of anhydrous sulfanilic acid to the solution and heat it until it dissolves. Cool the solution to room temperature, add 0.16 g of sodium nitrite, and stir the mixture until it is completely dissolved. Cool the solution in an ice-water bath for 5–10 minutes with stirring, or until the temperature is below 10°C. When the solution is completely cooled, add 0.5 mL of concentrated hydrochloric acid to the flask. The diazonium salt of sulfanilic acid should separate as a finely divided white precipitate after 1–2 minutes. Keep this suspension in an ice bath until it is to be used.

Methyl Orange. Combine 0.28 mL of *N,N*-dimethylaniline and 0.2 mL of glacial acetic acid in a test tube. Using a Pasteur pipette, add this solution to the cooled suspension of diazotized sulfanilic acid in the 50-mL Erlenmeyer flask. Stir the mixture vigorously. If the magnetic stirring bar does not move, you will need to use a glass stirring rod or spatula. Within a few minutes, a red precipitate of helianthin should form. Keep this mixture cooled for about 10 minutes in an ice bath. Next, add 3.0 mL of 10% aqueous sodium hydroxide. Do this slowly, with stirring, while you continue to cool the mixture in the ice bath. Using a glass rod, transfer a drop of the supernatant liquid to a strip of pH paper to determine if the solution is basic. If the mixture is too thick to perform this test, heat the mixture for a few

minutes first. If the solution is not basic, add more sodium hydroxide solution and check it again.

Heat the basic solution to boiling for 5–10 minutes to dissolve most of the newly formed methyl orange. It will not be possible to dissolve all of it. When most of the methyl orange has dissolved, add 1.0 g of sodium chloride and stir. Allow the mixture to cool to room temperature and then place in an ice bath to complete crystallization. Collect the product by vacuum filtration (Technique 8, Section 8.3) using a small Büchner funnel. Rinse the Erlenmeyer flask with two 2-mL ice-cold portions of saturated aqueous sodium chloride solution, and wash the filter cake of material you just collected with these rinse solutions.

To further purify the product, transfer the filter cake and paper to a 50-mL beaker containing about 12 mL of boiling water. Remove the paper with forceps. Maintain the mixture at a gentle boil for a few minutes, stirring it constantly with a glass stirring rod or microspatula. Not all of the dye will dissolve, but the salts with which it is contaminated will dissolve. Allow the solution to cool to room temperature. Cool the mixture in an ice bath. When it is cold, collect the product by vacuum filtration, using a small Büchner funnel. Allow the product to dry until the next laboratory period, weigh it, and calculate the percentage yield. Since salts do not generally have well-defined melting points, the melting-point determination should not be attempted. The methyl orange synthesized in this experiment can be used in Experiment 48, Formulation of a Paint and Art Project.

Tests (optional). Obtain a strip of a multifiber test fabric from your instructor. The cloth contains bands of different fabrics. Your instructor will identify which fabrics are found in your strip. Prepare a dye bath by dissolving 0.050 g of methyl orange (either your crude material or purified methyl orange can be used) in 30 mL of water to which 1 mL of a 15% aqueous sodium sulfate solution and 1 drop of concentrated sulfuric acid have been added. Heat the solution to just below its boiling point, immerse the fabric in the bath for 5–10 minutes. Remove the fabric, rinse it well, and note the results.

Make the dye bath basic by adding solid sodium carbonate. Then add a solution of sodium dithionite (sodium hydrosulfite) until the color of the bath is discharged. Add a slight excess of sodium dithionite. Now place the very end portion of the dyed fabric in the bath for a few minutes. Note the result.

Indicator Action (optional). Dissolve a few crystals of methyl orange in a small amount of water in a test tube. Alternately add a few drops of dilute hydrochloric acid and a few drops of dilute sodium hydroxide solution until the color changes are apparent in each case.

QUESTIONS

1. Why does the dimethylaniline couple with the diazonium salt at the *para* position of the ring?
2. Give a mechanism for producing a phenol from the diazonium salt that was prepared from sulfanilic acid.
3. What would be the result if copper (I) chloride were added to the diazonium salt prepared in the reaction?
4. The diazonium coupling reaction is an electrophilic aromatic substitution reaction. Give a mechanism that clearly indicates this fact.
5. Immediately after the coupling reaction in this experiment, a proton transfer occurs. The proton is transferred from the dimethylamino group to the azo linkage. Why is the latter protonated form lower in energy than the product formed initially?


Preparation of Indigo¹

Aldol condensation

Vat dyes

Ultraviolet spectroscopy

Indigo can be prepared by a series of reactions beginning with 2-nitrobenzaldehyde and acetone in the presence of base:

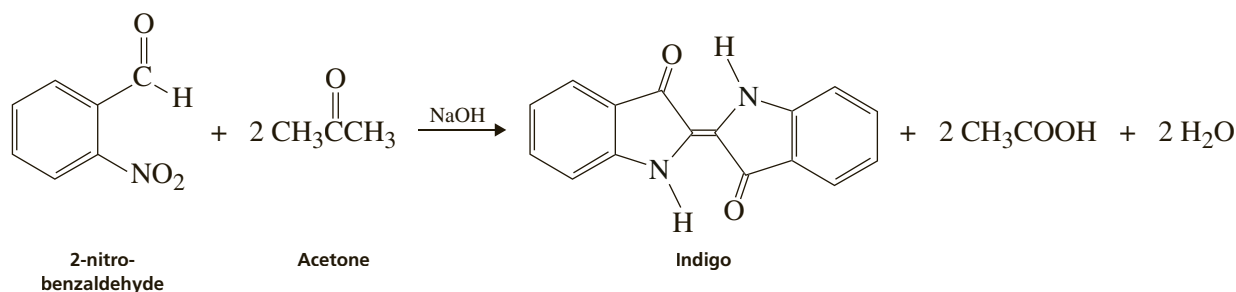

The first step in the reaction sequence is a crossed aldol condensation between 2-nitrobenzaldehyde and acetone. Subsequent steps in the formation of indigo are not fully understood.

After synthesizing indigo, you will dye a piece of multifiber fabric with your indigo. Since indigo is not water soluble, it must be first be reduced with sodium dithionite to form a water soluble, yellow form called leucoindigo. After the leucoindigo has been absorbed into the fibers of the cloth, the fabric is rinsed and allowed to dry. The reduced form of indigo then undergoes oxidation back to the insoluble form with the characteristic color of indigo. (See "Synthetic Dyes" essay, preceding Experiment 46.)

The ultraviolet spectrum of indigo will be obtained in this experiment. From this spectrum, it is possible to determine two defining characteristics about indigo: the wavelength maxima (λ_{\max}) and the molar absorptivity. The λ_{\max} is the wavelength at which the substance absorbs the most radiation (has the greatest absorption). The molar absorptivity ϵ is a wavelength-dependent physical property of the molecule undergoing electronic transition. It relates to the strength of the absorption at the specified wavelength λ_{\max} and will vary at other wavelengths. λ_{\max} is determined by looking at the spectrum and finding the wavelength where the absorption is greatest. The molar absorptivity is determined by using Beer's law:

$A = \epsilon cl$ for a given wavelength

A = absorbance (determined from the spectrum)

ϵ = molar absorptivity

c = molar concentration of solute

l = length of sample cell (cm)

¹ Parts of this experiment were adapted from "Ancient Dyes: Synthesis of Indigo," written by Amanda Murphy, Department of Chemistry, Western Washington University.

REQUIRED READING

Review: Technique 8 Filtration, Sections 8.1–8.6
New: Essay Synthetic Dyes

SPECIAL INSTRUCTIONS

Since indigo will likely stain the Büchner funnel permanently, it is best to reserve several funnels exclusively for this experiment. Indigo will stain your skin and clothes! You should wear gloves throughout the experiment. The synthesis of indigo requires ether to wash the product on the Büchner funnel. Ether is extremely flammable. Be certain that no open flames of any sort are in your vicinity when you are using ether.

NOTES TO THE INSTRUCTOR

A fresh bottle of sodium dithionite should be used. If students are performing the UV Characterization part of the experiment, the analytical balance should be accurate to 0.0001 g. For the dyeing part of the experiment, we use Multi-Fiber Ribbon (Item #SM-6) from Educational Innovation, 1-888-912-7434.

The indigo synthesized in this experiment can be used in Experiment 48, Formulation of a Paint and Art Project.

SUGGESTED WASTE DISPOSAL

Any organic waste should be discarded in an organic waste container designated for the disposal of *nonhalogenated* wastes.

PROCEDURE

Synthesis of Indigo. Add 0.5 g of *o*-nitrobenzaldehyde and 5.0 mL of acetone to a 50-mL beaker along with a magnetic stir bar. Stir the mixture until the solid dissolves. While continuing to stir, add 5.0 mL of distilled water to the beaker, which should result in a precipitate. With stirring, add 2.5 mL of 1.0 M sodium hydroxide dropwise. The solid will dissolve, and the solution should gradually darken as solid indigo is formed. Continue to stir for 10 minutes. Cool the mixture in an ice-bath for about 5 minutes and collect the indigo in a small Büchner funnel using vacuum filtration.² Wash the solid successively with 10 mL of water, 5 mL of ethanol, and 5 mL of ether. Allow the solid to dry on the Büchner funnel for 10 minutes and transfer the indigo to a tared watch glass. When the indigo has dried completely, determine the mass. Calculate the percent yield of indigo. The indigo synthesized in this experiment can also be used in Experiment 48, Formulation of a Paint and Art Project.

UV Characterization. Using an analytical balance, weigh approximately 1.0 mg of your indigo and place it in a 20-mL beaker. The weight need not be 1.0 mg exactly, but record the actual weight used to the nearest 0.1 mg. Add a magnetic stirring bar along with 10 mL of dimethylsulfoxide (DMSO) and stir until completely dissolved. Using an automatic pipette, transfer 0.100 mL of this solution

²Since indigo will permanently stain the Büchner funnel, your instructor may provide a funnel that is only used in this experiment.

to a cuvette. Dilute this with 1.90 mL of DMSO. Use a Pasteur pipette to mix the solution thoroughly by pulling up the liquid into the pipette and expelling several times. Record the UV spectrum, as directed by your instructor. Be sure to record the λ_{max} and the absorbance at this wavelength. Calculate the concentration of indigo in the final solution in moles per liter (you must take into account the dilution that you made). Using Beer's Law, calculate the molar absorptivity for your indigo sample (units in $M^{-1}\text{cm}^{-1}$). Compare your λ_{max} value and the molar absorptivity you calculated with the literature values.

Reduction of Indigo and Dyeing A Cloth. Obtain a strip of Multifiber Fabric from your instructor. This cloth contains acetate, SEF, Arnel, bleached cotton, Creslan 61, Dacron 54, Dacron 64, nylon 6.6, Orlon 75, spun silk, polypropylene, viscose, and wool in sequence. This fabric is used to demonstrate how the color of indigo appears on different fibers. Add 0.050 g of your dye and 0.50 g of sodium dithionite to a 50-mL filter flask. Using a glass stirring rod, grind the solids together. Add 5.0 mL of 1.0 M sodium hydroxide, 15 mL of water, and a magnetic stir bar to the flask. Place a rubber or cork stopper in the top, but leave the side-arm open. Heat the mixture on a hot plate with stirring. Initially the mixture will be a deep blue color, and you will not be able to see through the opaque liquid. Continue heating and stirring until the mixture has become a clear yellow. This color change indicates when the indigo has become reduced.³ When the liquid has become clear, add your fabric to the beaker. Allow the fabric to soak for 5 minutes and then remove it with a forceps. Immediately run the fabric under cold water until the color intensity on the fabric is constant and no longer decreases. This should take a minute or less. Then remove the fabric with your forceps and place the fabric on a paper towel for drying. Note the results in your notebook.

REFERENCE

- McKee, J. R.; Zanger, M. A Microscale Synthesis of Indigo: Vat Dyeing. *J. Chem. Edu.* **1991**, *68*, A242.
Boykin, D. W. A Convenient Apparatus for Small-Scale Dyeing with Indigo. *J. Chem. Edu.* **1998**, *75*, 769.

QUESTIONS

1. The first step in the reaction sequence for the formation of indigo involves a crossed aldol condensation between 2-nitrobenzaldehyde and acetone. Give the structure of the product of this first step.

³ If the liquid does not become clear or if you are unsure if it has become clear, add an additional 0.2–0.3 g of sodium dithionite and continue to heat the mixture for a few minutes.

Formulation of a Paint and Art Project¹

Formulation of a paint

Dyes

In this experiment, you will formulate a paint using either dyes that have been made in previous experiments (methyl orange in Experiment 46 and indigo in Experiment 47) or other dyes provided for this experiment. Using the paint that you make, you will create a painting!

The paint will generally consist of the following components:

Gum Arabic, a branched polysaccharide composed of galactose, rhamnose, arabinose, and glucuronic acid, functions as a **binder** that holds the dye in a fine suspension and fixes the paint to the paper or canvas.

Glycerol (glycerine) functions as a **plasticizer** that keeps the paint soft and easier to apply to the paper or canvas. It also helps the paint retain moisture and not dry out.

Sugar or honey functions as a **humectant** to help the paint retain moisture.

Dyes are used as **colorants**.

Water or ethanol can be added as a **paint thinner**.

Ethanol can provide better **wetting** of the paint and will shorten the drying time.

Starch or dextrin functions as a **filler** to thicken the paint and to make it smoother.

Alconox, or a similar laboratory **detergent**, also can be added to make the paint smoother.

Titanium dioxide and barium sulfate are **brighteners**, white compounds that can be added to give a brighter tint to the dye.

Calcium carbonate is used as a **thickener** to make the paint opaque.

Dyes are used in this experiment as **colorants** rather than the usual commercial pigments found in paints purchased at retail stores. Commercial paint pigments may be organic or inorganic compounds, but they are insoluble in the medium in which they are dispersed. The commercial paint is formulated to cover or “coat” the surface to which it is applied. The distinct particles of the pigment impart color and opacity to the paint. In contrast, dyes are organic compounds and are usually soluble in the medium (solvent) in which they are dispersed. Therefore, there are no visible particles of the dye, and the medium will retain its transparency. In normal use, dyes either react with the substrate or penetrate the fibers, such as those of wool or cotton, and impart color into the material being dyed. In this experiment calcium carbonate is used to make the paint more opaque, since the dyes we are

¹This experiment is based on a similar one developed by James Patterson, North Seattle College, Seattle, WA. The recipe for the paint is partly based on recipes developed by Patrician S. Hill, Professor of chemistry, Millersville University, Millersville, PA; and John Bordley, Professor emeritus of chemistry, Sewanee: The University of the South, Sewanee, TN.

using are fully transparent. If a transparent paint (more like a watercolor paint) is desired, the calcium carbonate can be omitted.

REQUIRED READING

New: Essay Synthetic Dyes

NOTES TO THE INSTRUCTOR

At North Seattle College and the University of Washington, we do this experiment on the last laboratory day of the quarter using dyes that have been prepared in previous experiments. Students enjoy creating a painting in their organic lab! Methyl orange (Experiment 46) and indigo (Experiment 47) can be used to color the paints. Other dyes can be provided to make the paintings more interesting. We select a panel of “experts” to select the best paintings. In almost every class, there are several paintings that are outstanding!

To make up the gum Arabic solution, place 56 g of powdered gum Arabic in a beaker and pour 118 mL of boiling water over the gum. Stir until it is dissolved. If lumps are still present, allow the mixture to stand for a while and stir again. Do not heat more after you have mixed the gum Arabic and boiling water.

Alternatively, a premade solution of gum Arabic can be purchased. We use a gum Arabic from Winsor & Newton (Water Color, Medium) which is readily available from most art supply stores.

Other dyes that we have made available to the students include Amaranth Red and Chicago Sky Blue.

We use watercolor paper for the paintings. The paper can be cut to different sizes with a paper cutter, depending on how large a painting each student wants to make.

WASTE DISPOSAL

Dispose of any leftover paint in the waste container for nonhalogenated organic waste.

PROCEDURE

The goal of this experiment is to formulate (make) a few milliliters of a paint that has a good color and consistency. You can use one of the dyes that you synthesized in a previous experiment (methyl orange or indigo), or you can use another dye that will be available in the lab. Using the paint that you make, you will create a painting. If you want to use more than one color in your painting, you can make up additional paints or use one created by one of your classmates.

Some of the components of paint listed in the introduction to this experiment are necessary and some are optional. This will be clearly explained in the procedure below. Formulating a paint can be done by following a precise recipe, but it can be more satisfying to experiment a little to see if you can come up with a paint that is most suitable for the kind of painting that you want to do. For example, you might want to use a paint that is more transparent, similar to a watercolor paint. In this case, you would not use any calcium carbonate or not very much. On the other hand, if you prefer a more opaque paint, you would use more calcium carbonate. Depending on what you are going to paint, you may want to formulate a paint that is thinner or thicker. To make a thinner paint, you would add more solvent (water or ethyl alcohol). To make a thicker paint, you would add more starch or dextrin.

Everyone should follow the instructions in this paragraph. Weigh 0.05–0.1 g of your dye and grind up finely in a mortar and pestle. Add 3–4 mL gum Arabic solution and grind this to get good mixing with the dye. Add 3–4 drops of honey and 3–4 drops of glycerine to the mortar. Mix thoroughly until you have a homogeneous mixture. If you desire a paint that has a more intense color, you should add more dye at this point. If you want to add another dye to achieve a new color, this is a good time to do this, too. Be sure to mix thoroughly after you have added more dye.

You can now adjust the appearance of your paint by adding the following components. Determine how much of each to add based on the consistency and appearance that you wish to achieve. In each case you should grind the mixture with the mortar and pestle to get good mixing. Alternatively, you can use your spatula to mix the ingredients.

1. Add solid calcium carbonate to make your paint more opaque.
2. Add water or ethanol to make the paint thinner. Ethanol will provide better “wetting” of the paint and will shorten the drying time.
3. Add starch or dextrin to thicken the paint and to make it smoother.
4. Add 1–2 granules of Alconox (or other laboratory detergent) to make the paint smoother.
5. Add solid titanium dioxide or barium sulfate to give a brighter tint to the paint.

When you have a paint that you are happy with, obtain a piece of watercolor paper of the desired size and a paintbrush. Paint your masterpiece! Your instructor will tell you where to place your painting to dry when you have finished it.

Place any leftover paint in a vial and label with your name. Alternatively, dispose of any leftover paint in the waste container for nonhalogenated organic waste.

REFERENCES

- Daines, T. L.; Morse, K. W. The Chemistry Involved in the Preparation of a Paint Pigment. *J. Chem. Educ.* **1976**, 53 (2) (Feb), 117–118.
- Mayer, R. *The Artist's Handbook*; Viking (the Penguin Group): New York, 1991.

INTERNET SITES

- <http://www.specialchem4coatings.com/tc/color-handbook/?id=physicalform>
<http://www.handprint.com/HP/WCL/pigmt1.html>

Local Anesthetics

Local anesthetics, or “painkillers,” are a well-studied class of compounds. Chemists have shown their ability to study the essential features of a naturally occurring drug and to improve on them by substituting totally new, synthetic surrogates. Often such substitutes are superior in desired medical effects and have fewer unwanted side effects or hazards.


The coca shrub (*Erythroxylon coca*) grows wild in Peru, specifically in the Andes Mountains, at elevations of 1,500 to 6,000 feet above sea level. The natives of South America have long chewed these leaves for their stimulant effects. Leaves of the coca shrub have even been found in pre-Inca Peruvian burial urns. Chewing the leaves brings about a definite sense of mental and physical well-being and the power to increase endurance. For chewing, the Indians smear the coca leaves with lime and roll them. The lime, Ca(OH)_2 , apparently releases the free alkaloid components; it is remarkable that the Indians learned this subtlety long ago by some empirical means. The pure alkaloid responsible for the properties of the coca leaves is **cocaine**.

The amounts of cocaine the Indians consume in this way are extremely small. Without such a crutch of central-nervous-system stimulation, the natives of the Andes would probably find it more difficult to perform the nearly Herculean tasks of their daily lives, such as carrying heavy loads over the rugged mountainous terrain. Unfortunately, overindulgence can lead to mental and physical deterioration and eventually an unpleasant death.

The pure alkaloid in large quantities is a common drug of addiction. Sigmund Freud first made a detailed study of cocaine in 1884. He was particularly impressed by the ability of the drug to stimulate the central nervous system, and he used it as a replacement drug to wean one of his addicted colleagues from morphine. This attempt was successful, but unhappily, the colleague became the world’s first known cocaine addict.

An extract from coca leaves was one of the original ingredients in Coca-Cola. However, early in the present century, government officials, with much legal difficulty, forced the manufacturer to omit coca from its beverage. The company has managed to this day to maintain the *coca* in its trademarked title, even though “Coke” contains none.

Our interest in cocaine lies in its anesthetic properties. The pure alkaloid was isolated in 1862 by Niemann, who noted that it had a bitter taste and produced a queer numbing sensation on the tongue, rendering it almost devoid of sensation. (Oh, those brave, but foolish chemists of yore who used to taste everything!) In 1880, Von Anrep found that the skin was made numb and insensitive to the prick of a pin when cocaine was injected subcutaneously. Freud and his assistant Karl Koller, having failed at attempts to rehabilitate morphine addicts, turned to a study of the anesthetizing properties of cocaine. Eye surgery is made difficult by involuntary reflex movements of the eye in response to even the slightest touch. Koller

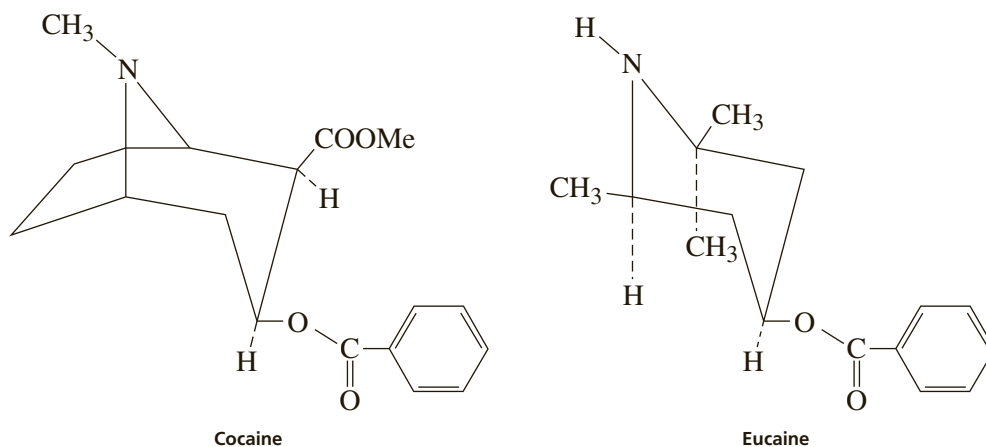


found that a few drops of a solution of cocaine would overcome this problem. Not only can cocaine serve as a local anesthetic but also it can be used to produce **mydriasis** (dilation of the pupil). The ability of cocaine to block signal conduction in nerves (particularly of pain) led to its rapid medical use in spite of its dangers. It soon found use as a “local” in both dentistry (1884) and in surgery (1885). In this type of application, it was injected directly into the particular nerves it was intended to deaden.

Soon after the structure of cocaine was established, chemists began to search for a substitute. Cocaine has several drawbacks for wide medical use as an anesthetic. In eye surgery, it also produces mydriasis. It can also become a drug of addiction. Finally, it has a dangerous effect on the central nervous system.


The first totally synthetic substitute was *eucaïne*. It was synthesized by Harries in 1918 and retains many of the essential skeletal features of the cocaine molecule. The development of this new anesthetic partly confirmed the portion of the cocaine structure essential for local anesthetic action. The advantage of eucaïne over cocaine is that it does not produce mydriasis and is not habit forming. Unfortunately, it is highly toxic.

A further attempt at simplification led to piperocaine. The molecular portion common to cocaine and eucaïne is outlined by dotted lines in the structure shown here. Piperocaine is only a third as toxic as cocaine itself.

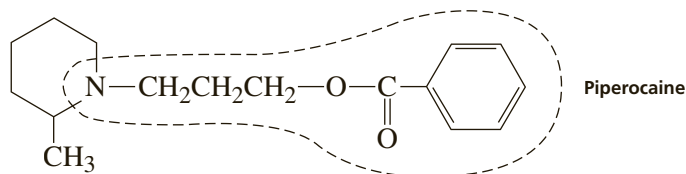


The most successful synthetic for many years was the drug procaine, known more commonly by its trade name Novocain (see table). Novocain is only a fourth as toxic as cocaine, giving a better margin of safety in its use. The toxic dose is almost 10 times the effective amount, and it is not a habit-forming drug.

Over the years, hundreds of new local anesthetics have been synthesized and tested. For one reason or another, most have not come into general use. The search for the perfect local anesthetic is still under way. All the drugs found to be active have certain structural features in common. At one end of the molecule is an aromatic ring. At the other is a secondary or tertiary amine. These two essential features are separated by a central chain of atoms usually one to four units long. The aromatic part is usually an ester of an aromatic acid. The ester group is important to the bodily detoxification of these compounds. The first step in deactivating them is a hydrolysis of this ester linkage, a process that occurs in the bloodstream.

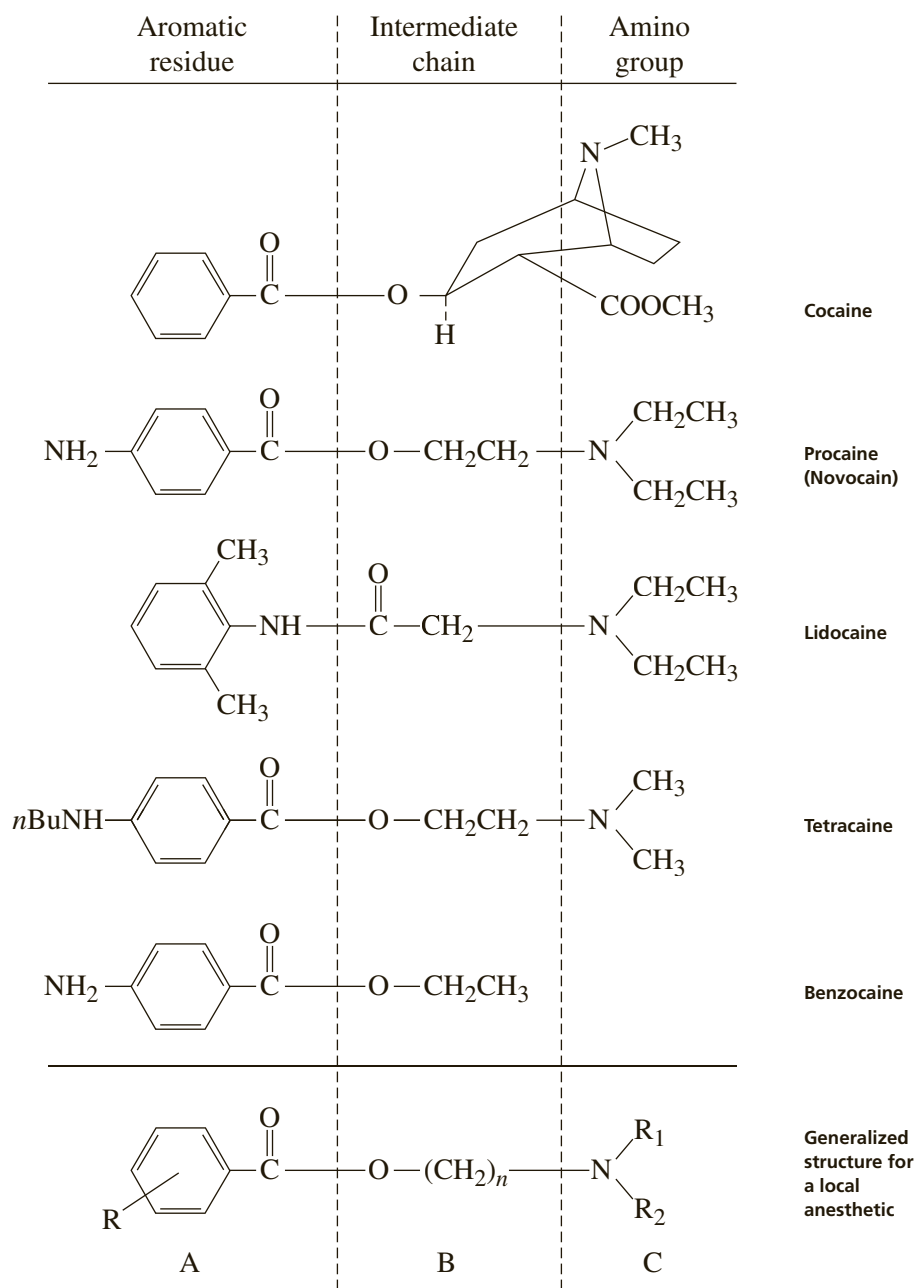

Local anesthetics.

Compounds that do not have the ester link are both longer lasting in their effects and generally more toxic. An exception is lidocaine, which is an amide. The tertiary amino group is apparently necessary to enhance the solubility of the compounds in the injection solvent. Most of these compounds are used in their hydrochloride salt forms, which can be dissolved in water for injection.

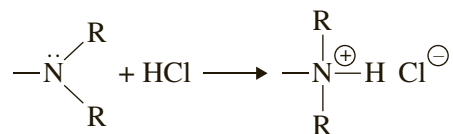

Benzocaine, in contrast, is active as a local anesthetic but is not used for injection. It does not suffuse well into tissue and is not water soluble. It is used primarily in skin preparations, in which it can be included in an ointment or salve for direct application. It is an ingredient of many sunburn-relief preparations.

How these drugs act to stop pain conduction is not well understood. Their main site of action is at the nerve membrane. They seem to compete with calcium at some receptor site, altering the permeability of the membrane and keeping the nerve slightly depolarized electrically.

REFERENCES


- Block, J. H.; Beale, J. M. *Wilson and Gisvold's Textbook of Organic Medicinal and Pharmaceutical Chemistry*, 11th ed.; Lippincott, Williams, and Wilkins: Philadelphia, 2003.
- Catterall, W. A.; Mackie, K. *Local Anesthetics*. In *Goodman and Gilman's, The Pharmacological Basis of Therapeutics*, 11th ed.; McGraw-Hill: New York, 2006.
- Hart, C.; Ksir, C., *Drugs, Society and Human Behavior*, 15th ed., McGraw-Hill, New York, 2012.
- Lenke, T. L.; Williams, D. A. *Foye's Principles of Medicinal Chemistry*, 6th ed.; Lippincott, Williams, and Wilkins: Philadelphia, 2008.
- Nagle, H.; Nagle, B. *Pharmacology: An Introduction*, 5th ed.; McGraw-Hill: Boston, 2005.
- Rang, H. P.; Dale, M. M.; Ritter J. M.; Flower, R.J. *Rang and Dale's Pharmacology*; Elsevier: China, 2007.
- Ray, O. S. "Stimulants and Depressants." In *Drugs, Society, and Human Behavior*, 3rd ed.; C. V. Mosby: St. Louis, 1983.
- Snyder, S. H. The Brain's Own Opiates. *Chem. Eng. News* 1977, (Nov 28), 26–35.
- Taylor, N. *Plant Drugs That Changed the World*; Dodd, Mead: New York, 1965; pp. 14–18.

Benzocaine

Esterification

Crystallization (mixed solvent method)

In this experiment, a procedure is given for the preparation of a local anesthetic, benzocaine, by the direct esterification of *p*-aminobenzoic acid with ethanol. At the instructor's option, you may test the prepared anesthetic on a frog's leg muscle.

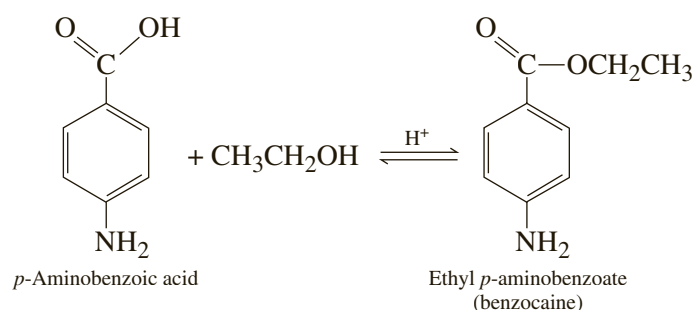

REQUIRED READING

- Review:* Technique 8 Filtration, Section 8.3
Technique 11 Crystallization, Sections 11.4 and 11.10
New: Essay Local Anesthetics

SPECIAL INSTRUCTIONS¹

Sulfuric acid is corrosive. Do not allow it to touch your skin.

SUGGESTED WASTE DISPOSAL

Dispose of all filtrates into the container designated for nonhalogenated organic solvents.

¹ *Note to the Instructor:* Benzocaine may be tested for its effect on a frog's leg muscle. See Instructor's Manual for instructions.

PROCEDURE

Running the Reaction

Place 0.120 g of *p*-aminobenzoic acid and 1.20 mL of absolute ethanol into a 3-mL conical vial. Add a magnetic spin vane and stir the mixture until the solid dissolves completely. While stirring, add 0.10 mL of concentrated sulfuric acid dropwise. A large amount of precipitate forms when you add the sulfuric acid, but this solid slowly dissolves during the reflux that follows. Attach a water-cooled condenser and heat the mixture at a gentle boil for 60–75 minutes with an aluminum block at about 105°C. Stir the mixture during this heating period.

Precipitation of Benzocaine

At the end of the reaction time, remove the apparatus from the aluminum block and allow the reaction mixture to cool for several minutes. Using a Pasteur pipette, transfer the contents of the vial to a small beaker containing 3.0 mL of water. When the liquid has cooled to room temperature, add a 10% sodium carbonate solution (about 1 mL needed) dropwise to neutralize the mixture. Stir the contents of the beaker with a stirring rod or spatula. After each addition of the sodium carbonate solution, extensive gas evolution (frothing) will be perceptible until the mixture is nearly neutralized. As the pH increases, a white precipitate of benzocaine is produced. When gas no longer evolves as you add a drop of sodium carbonate, check the pH of the solution and add further portions of sodium carbonate until the pH is about 8.


Collect the benzocaine by vacuum filtration, using a Hirsch funnel. Use three 1-mL portions of water to aid in the transfer and to wash the product in the funnel. Be sure that the solid is rinsed thoroughly with the water. After the product has dried overnight, weigh it, calculate the percentage yield, and determine its melting point. The melting point of pure benzocaine is 92°C.

Recrystallization and Characterization of Benzocaine


Although the product should be fairly pure, it may be recrystallized by the mixed solvent method using methanol and water (Technique 11, Section 11.10). Place the product in a Craig tube, add several drops of methanol, and, while heating the Craig tube in an aluminum block (60–70°C) and stirring the mixture with a microspatula, add methanol dropwise until all the solid dissolves. Add two to three additional drops of methanol and then add hot water dropwise until the mixture turns cloudy or a white precipitate forms. Add methanol again until the solid dissolves completely. Insert the inner plug of the Craig tube and allow the solution to cool slowly to room temperature. Complete the crystallization by cooling the mixture in an ice bath and collect the crystals by centrifugation (Technique 8, Section 8.8). Weigh the purified benzocaine and determine its melting point.

At the option of the instructor, obtain the infrared spectrum using the dry film method (Technique 25, Section 25.4) or by an alternative method (Technique 25, Section 25.5) and the NMR spectrum in CDCl_3 (Technique 26, Section 26.1)². Submit the sample in a labeled vial to the instructor.

² If a 60-MHz NMR spectrometer is used to determine the proton NMR of benzocaine, the amino protons may partially overlap the quartet in the ethyl group. If this is the case, a small amount of deuterated benzene can be added to shift the broad peak for the $-\text{NH}_2$ group away from the quartet: Carpenter, S. B.; R. H. Wallace. A Quick and Easy Simplification of Benzocaine's NMR Spectrum. *J. Chem. Educ.* **2006**, 83 (Apr), 637. A higher-field NMR spectrometer, such as obtained on a 300 MHz instrument, also avoids the overlap problem.

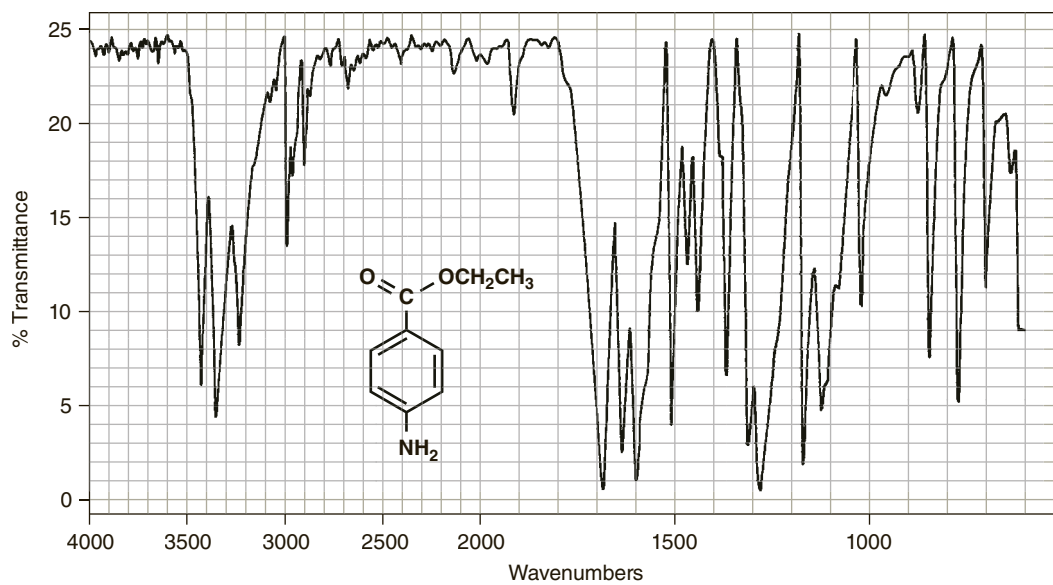

Infrared spectrum of benzocaine, KBr.

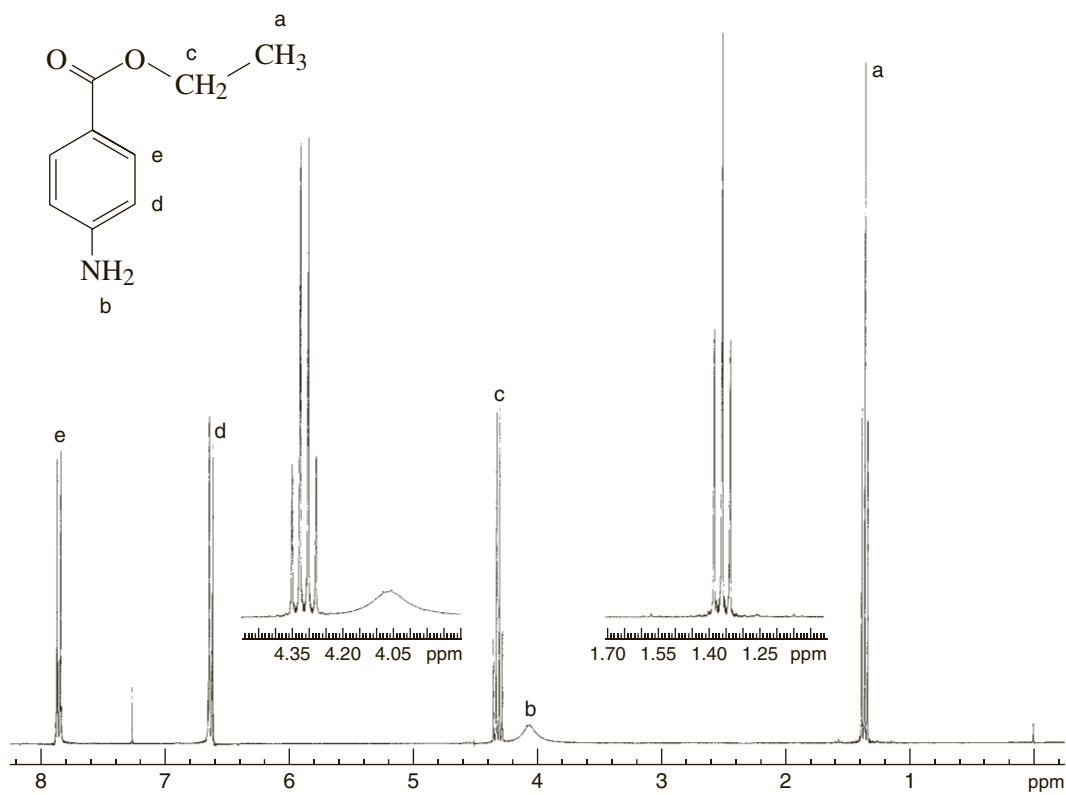300-MHz proton NMR spectrum of benzocaine (CDCl_3).

QUESTIONS

1. Interpret the infrared and NMR spectra of benzocaine.
2. What is the structure of the precipitate that forms after the sulfuric acid has been added?
3. When 10% sodium carbonate solution is added, a gas evolves. What is the gas? Give a balanced equation for this reaction.
4. Explain why benzocaine precipitates during the neutralization.
5. Refer to the structure of procaine in the table in the essay "Local Anesthetics." Using *p*-aminobenzoic acid, give equations showing how procaine and procaine monohydrochloride could be prepared. Which of the two possible amino functional groups in procaine will be protonated first? Defend your choice. (*Hint*: Consider resonance.)

Pheromones: Insect Attractants and Repellents

It is difficult for humans, who are accustomed to heavy reliance on visual and verbal forms of communication, to imagine that some forms of life depend primarily on the release and perception of *odors* to communicate with one another. Among insects, however, this is perhaps the chief form of communication. Many species of insects have developed a virtual “language” based on the exchange of odors. These insects have well-developed scent glands, often of several different types, which have as their sole purpose the synthesis and release of chemical substances. When these chemical substances, known as **pheromones**, are secreted by insects and detected by other members of the same species, they induce a specific and characteristic response.

TYPES OF PHEROMONES

Releaser pheromones: This type of pheromone produces an immediate behavioral response, but is quickly dissipated. Releaser molecules can attract mates from considerable distances, but the effect is short-lived.

Primer pheromones: Primer pheromones trigger a series of physiological changes in the recipient. In contrast to a releaser pheromone, a primer pheromone has a slower onset and a longer duration.

Recruiting or aggregation pheromones: This type of pheromone can attract individuals of both sexes of the same species.

Recognition pheromones: This type of pheromone allows members of the same species to recognize one another. This type of pheromone serves a similar function to recruiting pheromones.

Alarm pheromones: This type of substance is released when attacked by a predator. It can alert others to escape, or it can cause an aggressive response to members of the same species.

Territorial pheromones: These pheromones mark the boundaries of an organism’s territory. In dogs, these pheromones are present in the urine. Dogs can thus mark out their territory.

Trail pheromones: Ants deposit a trail of pheromones as they return to the nest from their source of food. This trail attracts other ants and serves as a guide to the food source. The pheromone must be continually renewed because the low-molecular-weight compounds evaporate rapidly.

Sex pheromones: Sex pheromones indicate the availability of the female for breeding purposes. Male animals also emit pheromones that convey information about their species. No confusion results!

It should be mentioned that there is some overlap of function in pheromones. The pheromones can assume multiple responses even though they may be categorized separately.

SEX ATTRACTANTS

Among the most important types of releaser pheromones are the sex attractants. **Sex attractants** are pheromones secreted by either the female or, less commonly, the male of the species to attract the opposite sex for the purpose of mating. In large concentrations, sex pheromones also induce a physiological response in the recipient (for example, the changes necessary to the mating act). Thus, they have a primer effect and so are misnamed.

Anyone who has owned a female cat or dog knows that sex pheromones are not limited to insects. Female cats or dogs widely advertise, by odor, their sexual availability when they are “in heat.” This type of pheromone is not uncommon to mammals. Some persons even believe that human pheromones are responsible for attracting certain sensitive males and females to one another. This idea is, of course, responsible for many of the perfumes now widely available. Whether or not the idea is correct cannot yet be established, but there are proven sexual differences in the ability of humans to smell certain substances. For instance, Exaltolide, a synthetic lactone of 14-hydroxytetradecanoic acid, can be perceived only by females or males after they have been injected with an estrogen. Exaltolide is very similar in overall structure to civetone (civet cat) and muskone (musk deer), which are two naturally occurring compounds believed to be mammalian sex pheromones.

Whether or not humans use pheromones as a means of attracting the opposite sex has never been completely established, although it is an active area of research. Humans, like other animals, emit odors from many parts of their bodies. Body odor consists of secretions from several types of skin glands, most of which are concentrated in the underarm region of the body. Do these secretions contain substances that might act as human pheromones?

Research has shown that a mother can correctly identify the odor of her newborn infant or older child by smelling clothing worn previously by the child and can distinguish the clothing from that worn by another child of the same age. Studies conducted over 30 years ago showed that the menstrual cycles of women who are roommates or close friends tend to converge over time. These and other similar investigations suggest that some forms of pheromone-like communication are possible in humans.

Recent studies have clearly identified a specialized structure, called the **vomer nasal organ**, in the nose. This organ appears to respond to a variety of chemical stimuli. In a recent article, researchers at the University of Chicago reported that when they wiped human body-odor secretions from one group of women under the noses of other women, the second group showed changes in their menstrual cycles. The cycles grew either longer or shorter, depending on where the donors were in their own menstrual cycles. The affected women claimed that they did not smell anything except the alcohol on the cotton pads. Alcohol alone had no effect on the women’s menstrual cycles. The timing of ovulation for the female test subjects was affected in a similar manner. Although the nature

of substances responsible for these effects has not yet been identified, clearly the potential for chemical communication regulating sexual function has been established in humans.

This effect has been described as the McClintock effect, named after the primary investigator, Martha McClintock, at the University of Chicago (see references: McClintock and Stern, 1971 and 1998). The McClintock effect, however, is still not firmly established and more recent studies and reviews of the McClintock research have called into question the result of the study (see references: Yang and Schank, 2006).

One of the first identified insect attractants belongs to the gypsy moth, *Lymantria dispar*. This moth is a common agricultural pest, and it was hoped that the sex attractant that females emit could be used to lure and trap males. Such a method of insect control would be preferable to inundating large areas with insecticides and would be species-specific. Nearly 50 years of work were expended in identifying the chemical substance responsible for the attractant's power. Early in this period, researchers found that an extract from the tail sections of female gypsy moths would attract males, even from a great distance. Experiments with the isolated gypsy moth pheromone demonstrated that the male gypsy moth has an almost unbelievable ability to detect extremely small amounts of the substance. He can detect it in concentrations lower than a few hundred *molecules* per cubic centimeter (about 10^{-19} – 10^{-20} g/cc)! When a male moth encounters a small concentration of pheromone, he immediately heads into the wind and flies upward in search of higher concentrations and the female. In only a mild breeze, a continuously emitting female can activate a space 300 ft high, 700 ft wide, and almost 14,000 ft (nearly 3 miles) long!

In subsequent work, researchers isolated 20 mg of a pure chemical substance from solvent extracts of the two extreme tail segments collected from each of 500,000 female gypsy moths (about 0.1 μ g/moth). This emphasizes that pheromones are effective in very minute amounts and that chemists must work with very small amounts to isolate them and prove their structures. It is not unusual to process thousands of insects to get even a minute sample of these substances. Extremely sophisticated analytical and instrumental methods, such as spectroscopy, must be used to determine the structure of a pheromone.

In spite of these techniques, the original researchers assigned an incorrect structure to the gypsy moth pheromone and proposed for it the name **gyplure**. Because of its great promise as a method of insect control, gyplure was soon synthesized. The synthetic material turned out to be totally inactive. After some controversy about why the synthetic material was incapable of luring male gypsy moths (see the References for the complete story), it was finally shown that the proposed structure for the pheromone (that is, the gyplure structure) was incorrect. The actual pheromone was found to be *cis*-7,8-epoxy-2-methyloctadecane, also named (7R,8S)-epoxy-2-methyloctadecane. This material was soon synthesized, found to be active, and given the name **disparlure**. In recent years, disparlure traps have been found to be a convenient and economical method for controlling the gypsy moth.

A similar story of mistaken identity can be related for the structure of the pheromone of the pink bollworm, *Pectinophora gossypiella*. The originally proposed structure was called **propylure**. Synthetic propylure turned out to be inactive. Subsequently, the pheromone was shown to be a mixture of two isomers of 7,11-hexadecadien-1-yl acetate, the *cis,cis* (7Z,11Z) isomer and the *cis,trans* (7Z,11E) isomer. It turned out to be quite easy to synthesize a 1:1 mixture of these two isomers, and the 1:1 mixture was named **gossylure**. Curiously, adding as little as 10% of either

of the other two possible isomers, *trans,cis* (7E,11Z) or *trans,trans* (7E,11E), to the 1:1 mixture greatly diminishes its activity, apparently masking it. Geometric isomerism can be important! The details of the gossyplure story can also be found in the References.

Both these stories have been partly repeated here to point out the difficulties of research on pheromones. The usual method is to propose a structure determined by work on *very tiny* amounts of the natural material. The margin for error is great. Such proposals are usually not considered “proved” until synthetic material is shown to be as biologically effective as the natural pheromone.

OTHER PHEROMONES

The most important example of a primer pheromone is found in honeybees. A bee colony consists of one queen bee, several hundred male drones, and thousands of worker bees, or undeveloped females. It has recently been found that the queen, the only female that has achieved full development and reproductive capacity, secretes a primer pheromone called the **queen substance**. The worker females, while tending the queen bee, continuously ingest quantities of the queen substance. This pheromone, which is a mixture of compounds, prevents the workers from rearing any competitive queens and prevents the development of ovaries in all other females in the hive. The substance is also active as a sex attractant; it attracts drones to the queen during her “nuptial flight.” The major component of the queen substance is shown in the figure.

Honeybees also produce several other important types of pheromones. It has long been known that bees will swarm after an intruder. It has also been known that isopentyl acetate induces a similar behavior in bees. Isopentyl acetate (Experiment 14) is an **alarm pheromone**. When an angry worker bee stings an intruder, she discharges, along with the sting venom, a mixture of pheromones that incites the other bees to swarm on and attack the intruder. Isopentyl acetate is an important component of the alarm pheromone mixture. Alarm pheromones have also been identified in many other insects. In insects less aggressive than bees or ants, the alarm pheromone may take the form of a **repellent**, which induces the insects to go into hiding or leave the immediate vicinity.

Honeybees also release **recruiting or trail pheromones**. These pheromones attract others to a source of food. Honeybees secrete recruiting pheromones when they locate flowers in which large amounts of sugar syrup are available. Although the recruiting pheromone is a complex mixture, both geraniol and citral have been identified as components. In a similar fashion, when ants locate a source of food, they drag their tails along the ground on their way back to the nest, continuously secreting a trail pheromone. Other ants follow the trail to the source of food.

In some species of insects, **recognition pheromones** have been identified. In carpenter ants, a caste-specific secretion has been found in the mandibular glands of the males of five different species. These secretions have several functions, one of which is to allow members of the same species to recognize one another. Insects not having the correct recognition odor are immediately attacked and expelled from the nest. In one species of carpenter ant, methyl anthranilate has been shown to be an important component of the recognition pheromone.

We do not yet know all the types of pheromones that any given species of insect may use, but it seems that as few as 10 or 12 pheromones could constitute a “language” that could adequately regulate the entire life cycle of a colony of social insects.

INSECT REPELLENTS


Currently, the most widely used **insect repellent** is the synthetic substance *N,N*-diethyl-*m*-toluamide (see Experiment 50), also called Deet. It is effective against fleas, mosquitoes, chiggers, ticks, deerflies, sandflies, and biting gnats. A specific repellent is known for each of these types of insects, but none has the wide spectrum of activity that this repellent has. Exactly why these substances repel insects is not yet fully understood. The most extensive investigations have been carried out on the mosquito.












Originally, many investigators thought that repellents might simply be compounds that provided unpleasant or distasteful odors to a wide variety of insects. Others thought that they might be alarm pheromones for the species affected, or that they might be the alarm pheromones of a hostile species. Early research with the mosquito indicates that at least for several varieties of mosquitoes, none of these is the correct answer.

Mosquitoes seem to have hairs on their antennae that are receptors enabling them to find a warm-blooded host. These receptors detect the convection currents arising from a warm and moist living animal. When a mosquito encounters a warm and moist convection current, the mosquito moves steadily forward. If it passes out of the current into dry air, it turns until it finds the current again. Eventually, it finds the host and lands. Repellents cause a mosquito to turn in flight and become confused. Even if it should land, it becomes confused and flies away again.

Researchers have found that the repellent prevents the moisture receptors of the mosquito from responding normally to the raised humidity of the subject. At least two sensors are involved, one responsive to carbon dioxide and the other responsive to water vapor. The carbon dioxide sensor is activated by the repellent, but if exposure to the chemical continues, adaptation occurs, and the sensor returns to its usual low output of signal. The moisture sensor, on the other hand, simply seems to be deadened, or turned off, by the repellent. Therefore, mosquitoes have great difficulty in finding and interpreting a host when they are in an environment saturated with repellent. They fly right through warm and humid convection currents as if the currents did not exist. Only time will tell if other biting insects respond likewise.

Until now, the mechanism of action of insect repellents on molecular targets remained unknown. However, Leslie Vooshall and colleagues at Rockefeller University reported in the March 2008 issue of *Science* that they had identified the molecular targets for the repellent, *N,N*-diethyl-*meta*-toluamide (DEET). They reported that DEET inhibits mosquito and fruit fly olfactory receptors that form a complex with a required olfactory co-receptor, OR83b. In effect, DEET inhibits behavioral attraction by masking the host odor in humans. Now that it is known how DEET affects receptors, new insect repellants may be developed that are safer and more effective, especially for young children.

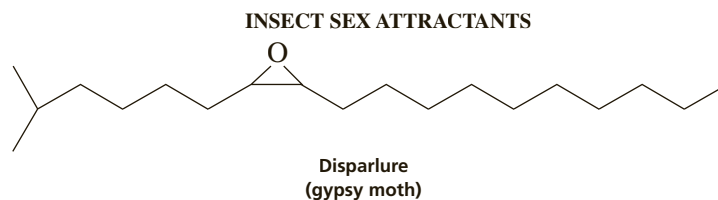

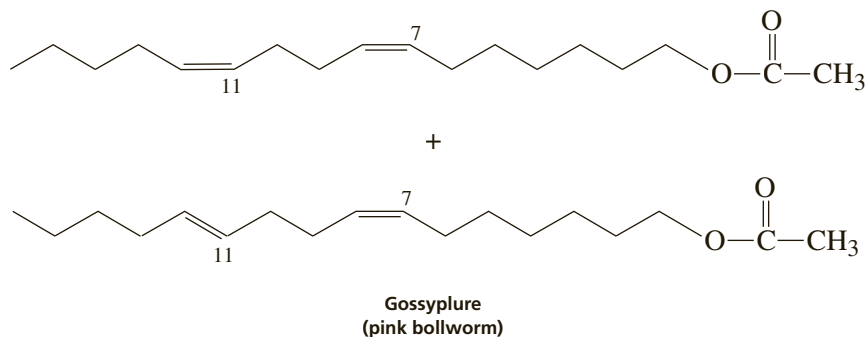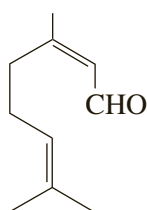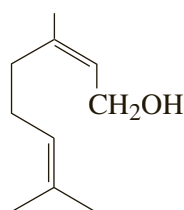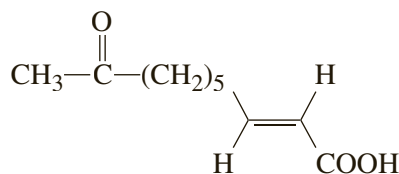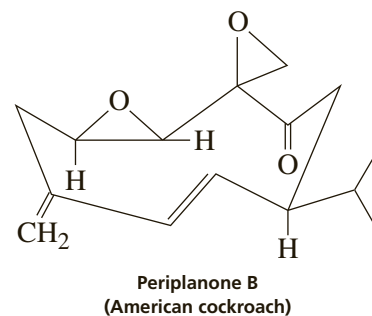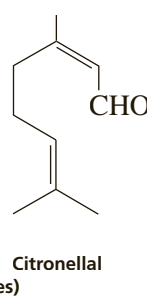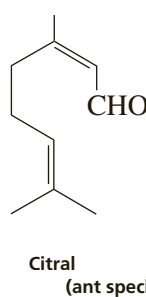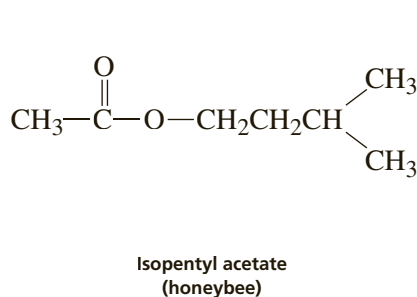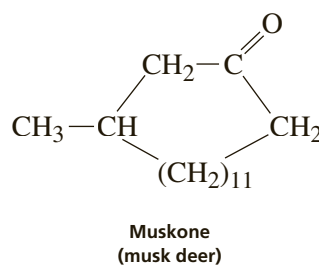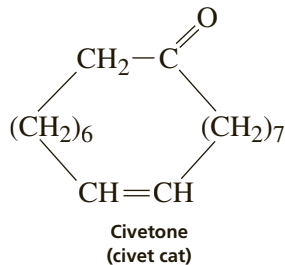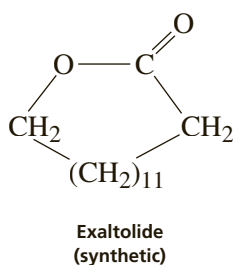
RECRUITING PHEROMONES

PRIMER PHEROMONE

ALARM PHEROMONES

MAMMALIAN PHEROMONES (?)

REFERENCES

- Agosta, W. C. Using Chemicals to Communicate. *J. Chem. Educ.* **1994**, 71 (Mar), 242.
- Batra, S. W. T. Polyester-Making Bees and Other Innovative Insect Chemists. *J. Chem. Educ.* **1985**, 62 (Feb), 121.
- Ditzen, M.; Pellegrino, M.; Vosshall, L. B. Insect Odorant Receptors Are Molecular Targets of the Insect Repellent DEET. *Science* **2008**, 319, 1838–42.
- Katzenellenbogen, J. A. Insect Pheromone Synthesis: New Methodology. *Science* **1976**, 194 (Oct 8), 139.

- Keeling, C. I.; Slessor, K. N.; Higo, H. A.; Winston, M. L. New Components of the Honey Bee (*Apis mellifera* L.) Queen Retinue Pheromone. *Proc. Natl. Acad. Sci.* **2003**, *100* (8) (Apr 15), 4486.
- Kohl, J. V.; Atzmueller, M.; Fink, B.; Grammer, K. Human Pheromones: Integrating Neuroendocrinology and Ethology. *Neuroendocrinol. Lett.* **2001**, *22* (5), 319–31.
- Leonhardt, B. A. Pheromones. *CHEMTECH* **1985**, *15* (Jun), 368.
- Liberles, S. D.; Buck, L. B. A Second Class of Chemosensory Receptors in the Olfactory Epithelium. *Nature* **2006**, *442*, 645–50.
- Prestwick, G. D. The Chemical Defenses of Termites. *Sci. Am.* **1983**, *249* (Aug), 78.
- Silverstein, R. M. Pheromones: Background and Potential Use for Insect Control. *Science* **1981**, *213* (Sept 18), 1326.
- Stine, W. R. Pheromones: Chemical Communication by Insects. *J. Chem. Educ.* **1986**, *63* (Jul), 603.
- Villemin, D. Olefin Oxidation: A Synthesis of Queen Bee Pheromone. *Chem. Ind.* **1986** (Jan 20), 69.
- Wilson, E. O. Pheromones. **1963**, *208* (May), 100.
- Wilson, E. O.; Bossert, W. H. Chemical Communication Among Animals. *Recent Progr. Horm. Res.* **1963**, *19*, 673–716.
- Winston, M. L.; Slessor, K. N. The Essence of Royalty: Honey Bee Queen Pheromone. *Am. Sci.* **1992**, *80* (Jul–Aug), 374.
- Wood, W. F. Chemical Ecology: Chemical Communication in Nature. *J. Chem. Educ.* **1983**, *60* (Jul), 531.
- Wright, R. H. Why Mosquito Repellents Repel. *Sci. Am.* **1975**, *233* (Jul), 105.
- Wyatt, Tristram D. *Pheromones and Animal Behaviour: Communication by Smell and Taste*; Cambridge University Press: Cambridge, 2003.
- Yu, H.; Becker, H.; Mangold, H. K. Preparation of Some Pheromone Bouquets. *Chem. Ind.* **1989**, (Jan 16), 39.

Gypsy Moth

- Beroza, M.; Knipling, E. F. Gypsy Moth Control with the Sex Attractant Pheromone. *Science* **1972**, *177*, 19.
- Bierl, B. A.; Beroza, M.; Collier, C. W. Potent Sex Attractant of the Gypsy Moth: Its Isolation, Identification, and Synthesis. *Science* **1970**, *170* (3953), 87.

Pink Bollworm

- Anderson, R. J.; Henrick, C. A. Preparation of the Pink Bollworm Sex Pheromone Mixture, Gossyplure. *J. Am. Chem. Soc.* **1975**, *97* (15), 4327.
- Hummel, H. E.; Gaston, L. K.; Shorey, H. H.; Kaae, R. S.; Byrne, K. J.; Silverstein, R. M. Clarification of the Chemical Status of the Pink Bollworm Sex Pheromone. *Science* **1973**, *181* (4102), 873.

Cockroach

- Adams, M. A.; Nakanishi, K.; Still, W. C.; Arnold, E. V.; Clardy, J.; Persoon, C. J. Sex Pheromone of the American Cockroach: Absolute Configuration of Periplanone-B. *J. Am. Chem. Soc.* **1979**, *101*, 2495.
- Nojima, S.; Schall, C.; Webster, F. X.; Santangelo, R. G.; Roelofs, W. L. Identification of the Sex Pheromone of the German Cockroach. *Blatella germanica* *Science* **2005**, *307*, 1104.
- Still, W. C. (±)-Periplanone-B: Total Synthesis and Structure of the Sex Excitant Pheromone of the American Cockroach. *J. Am. Chem. Soc.* **1979**, *101*, 2493.
- Stinson, S. C. Scientists Synthesize Roach Sex Excitant. *Chem. Eng. News* **1979**, *57* (Apr 30), 24.

Spider

- Schulz, S.; Toft, S. Identification of a Sex Pheromone from a Spider. *Science* **1993**, *260* (Jun 11), 1635.

Silkworm

- Emsley, J. Sex and the Discerning Silkworm. *Foodweek* **1992**, *135* (Jul 11), 18.

Aphids

- Coghlan, A. Aphids Fall for Siren Scent of Pheromones. *Foodweek* **1990**, *127* (Jul 21), 32.

Snakes

Mason, R. T.; Fales, H. M.; Jones, T. H.; Pannell, L. K.; Chinn, J. W.; Crews, D. Sex Pheromones in Snakes. *Science* **1989**, 245 (Jul 21), 290.

Oriental Fruit Moth

Mithran, S.; Mamdapur, V. R. A Facile Synthesis of the Oriental Fruit Moth Sex Pheromone. *Chem. Ind.* **1986**, (Oct 20), 711.

Human

McClintock, M. K. Menstrual Synchrony and Suppression. *Nature* **1971**, 229, 244–45.

Stern, K.; McClintock, M. K. Regulation of Ovulation by Human Pheromones. *Nature* **1998**, 392, 177–79.

Weller, A. Communication through Body Odour. *Nature* **1998**, 392 (Mar 12), 126.

Yang, Z. J. C. Women Do Not Synchronize Their Menstrual Cycles. *Hum. Nat.* **2006**, 17 (4), 434–47.

Worms (postpones aging)

Butcher, R. A.; Fujita, M.; Schroeder, F. C.; Clardy, J. Small-Molecule Pheromones that Control Dauer Development in *Caenorhabditis Elegans*. *Nature Chem. Bio.* **2007**, 3, 420.

Rats (Anxiety)

Sniffing out Anxiety in Rattled Rats: Hexanal and 4-Methylpentanal. *Chemical and Engineering News* **2014**, 92 (Dec 22), 12, referencing an article in *Proc. Natl. Acad. Sci. USA* **2014**, DOI: 10.1073/pnas.1414710112.

INTERNET SITES

Wikipedia, <http://en.wikipedia.org/wiki/Pheromone>. This site describes the various types of pheromones.

Pherobase, the database of insect pheromones
<http://www.pherobase.com/>


The Pherobase database is an extensive compilation of behavior-modifying compounds listed in the various pheromone categories: aggregation, alarm, releaser, primer, territorial, trail, sex pheromones, and others. The database contains over 30,000 entries. Jmol images of molecules are shown. The molecules can be projected as either space-filling or wire-frame models. They can be rotated in three-dimensional space. In addition, the database includes mass spectral, NMR, and synthesis data for more than 2,500 compounds. This is a fun site!

N,N-Diethyl-*m*-toluamide: The Insect Repellent "OFF"


Preparation of an amide

Extraction

In this experiment, you will synthesize the active ingredient of the insect repellent "OFF," *N,N*-diethyl-*m*-toluamide. This substance belongs to the class of compounds called **amides**. Amides have the generalized structure

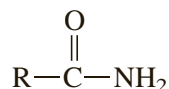

The amide to be prepared in this experiment is a disubstituted amide. That is, two of the hydrogens on the amide— NH_2 group have been replaced with ethyl groups. Amides cannot be prepared directly by mixing a carboxylic acid with an amine. If an acid and an amine are mixed, an acid–base reaction occurs, giving the conjugate base of the acid, which will not react further while in solution:

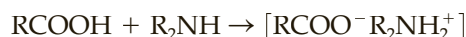


However, if the amine salt is isolated as a crystalline solid and strongly heated, the amide can be prepared:




Because of the high temperature required for this reaction, this is not a convenient laboratory method.

Amides are usually prepared via the acid chloride, as in this experiment. In step 1, *m*-toluic acid is converted to its acid chloride derivative using thionyl chloride (SOCl_2).

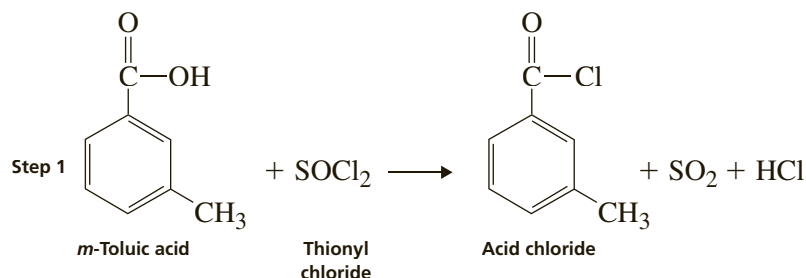

The acid chloride is not isolated or purified, and it is allowed to react directly with diethylamine in step 2. An excess of diethylamine is used in this experiment to react with the hydrogen chloride produced in step 2.

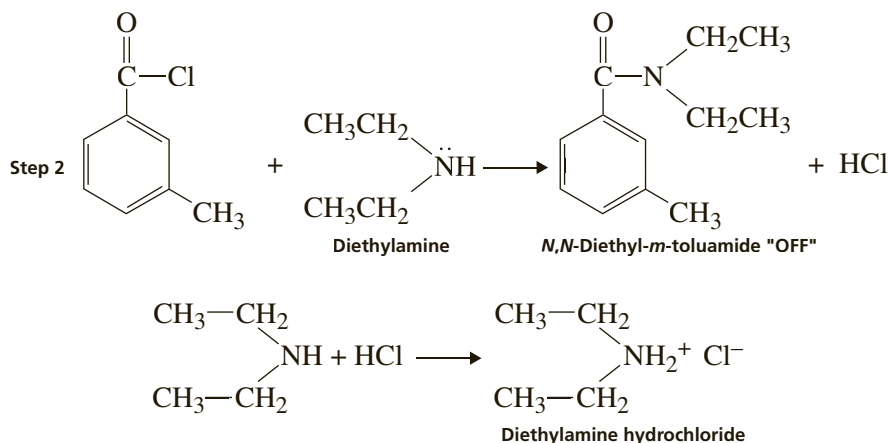

REQUIRED READING

- Review:* Technique 7 Reaction Methods, Sections 7.3 and 7.10
 Technique 12 Extractions, Separations, and Drying Agents,
 Sections 12.4, 12.8, 12.9, 12.11
- New:* Essay Pheromones: Insect Attractants and Repellents

SPECIAL INSTRUCTIONS

All equipment used in this experiment should be dry because thionyl chloride reacts with water to liberate HCl and SO₂. Likewise, *anhydrous* ether should be used because water reacts with both thionyl chloride and the intermediate acid chloride.

Thionyl chloride is a noxious and corrosive chemical and should be handled with care. If it is spilled on the skin, serious burns will result. Thionyl chloride and diethylamine must be dispensed *in the hood* from bottles that should be kept tightly closed when not in use. Diethylamine is also noxious and corrosive. In addition, it is quite volatile (BP 56°C) and must be cooled in a hood prior to use.

SUGGESTED WASTE DISPOSAL

All aqueous extracts should be poured into the waste bottle designated for aqueous waste.

PROCEDURE

Preparation of the Acid Chloride

Place 1.81 g (0.0133 mol) of *m*-toluic acid (3-methybenzoic acid, MW = 136.1) into a dry 25-mL round-bottom flask. Add 1 mL of anhydrous diethyl ether to wet the solid (it will not dissolve), and place a stir bar in the flask. In a hood, carefully add 2.0 mL of thionyl chloride (0.0275 mol, density = 1.64 g/mL, MW 118.9) from a plastic Pasteur pipette. Thionyl chloride is a nasty substance, so be careful not to breathe in the vapors! Add 5 drops of pyridine. At this point, you should observe a rapid reaction with evolution of gases. *Lightly stopper the flask.* The reaction will liberate sulfur dioxide and hydrogen chloride, so make sure that the flask is kept in a well-ventilated hood. Stir the mixture for about 10 minutes. During the course of the reaction period, the solid

m-toluic acid will slowly dissolve (react) with the thionyl chloride. Continue stirring until the solid has dissolved.

CAUTION

The thionyl chloride is kept in a hood. Do not breathe the vapors of this noxious and corrosive chemical. Use dry equipment when handling this material because it reacts violently with water. Do not get it on your skin.

Insert a piece of glass tubing through the rubber piece on the thermometer adapter and insert it into the neck of the 25-mL round-bottom flask. Remove the excess thionyl chloride under vacuum using an aspirator (with water trap!) or with the house vacuum system. The best way to remove the excess thionyl chloride is to swirl the flask, rather than using the magnetic stirring unit. Do not heat the mixture. The mixture will show obvious signs on boiling under the vacuum. You should see boiling action around the stir bar, accompanied by a little frothing. Continue to pull a vacuum on the flask until the boiling action ceases or nearly ceases. At that point, the volume should have been reduced. It may take about ½ hour to remove the excess thionyl chloride. Swirl the mixture continuously during this time to aid the evaporation process.

Preparation of the Amide

Prepare a solution of diethylamine in aqueous sodium hydroxide solution by adding 4 mL of diethylamine (0.430 mol, density = 0.71 g/mL, MW 73.1) from a plastic disposable Pasteur pipette to 15 mL of 10% aqueous sodium hydroxide solution in a 50-mL Erlenmeyer flask. Cool the mixture to 0°C in an ice-water bath. *Slowly* add the acid chloride mixture with a plastic Pasteur pipette to the cooled diethylamine/sodium hydroxide mixture, with swirling of the flask. *The reaction is violent, and a lot of smoke is observed.* Add the acid chloride in small portions over about a 5-minute period. Following the addition, swirl the mixture in the flask occasionally over a 10-min period to complete the reaction.

Isolation of the Amide

Pour the mixture into a separatory funnel using portions of 20 mL of diethyl ether to aid the transfer. Add the rest of the diethyl ether and shake the separatory funnel to extract the product from the aqueous mixture. Remove the lower aqueous layer and save it. Pour the ether layer out of the top of the funnel into an Erlenmeyer flask for temporary storage. Return the aqueous layer back into the separatory funnel and extract it with a fresh 20-mL portion of ether. Remove the aqueous layer and discard it. Pour the ether layer from the top of the separatory funnel and into the flask containing the first ether extract. Return the combined ether layers into the separatory funnel, and shake it with a 20-mL portion of saturated aqueous NaCl solution to do a preliminary drying of the ether layer. Remove the lower aqueous layer and discard it. Pour the ether solution from the top of the separatory funnel into a *dry* Erlenmeyer flask. Dry the ether layer with anhydrous magnesium sulfate. Decant the solution away from the drying agent through a piece of fluted filter paper into a *preweighed* 100-mL round-bottom flask. Remove the ether on a rotary evaporator or remove the ether under vacuum (see Technique 7, Figure 7.18C).¹


¹ You may need to obtain a 100-mL round-bottom flask from your instructor if using a rotary evaporator. Alternatively, if evaporating on the solvent using Technique 7, Figure 7.18C, then use a preweighed filter flask.

Reweigh the flask to determine the yield of the reddish-brown product. Yields are generally reasonable and exceed 80%.

Analysis of the Product

Determine the infrared spectrum of your product. The spectrum can be compared to the one reproduced in Figure 1. You may see a small amount of unreacted diethylamine appearing near 3400 cm^{-1} in your spectrum, with can be ignored.

At the option of the instructor determine the ^1H (proton) NMR spectrum of your product. The 500 MHz spectrum determined at 20°C shows interesting pattern for the ethyl groups attached to a nitrogen (Figure 2). The two methylene carbon atoms in the ethyl groups appear as a pair of broad peaks between 3.2 and 3.6 ppm, indicating non-equivalence. Notice that the peaks are broad and do not show up as quartets. Likewise, the two methyl carbon atoms in the ethyl groups appear as a pair of broad peaks between 1.0 and 1.3 ppm and do not show up as triplets. There is restricted rotation in amides resulting from resonance, leading to non-equivalence of the two ethyl groups:

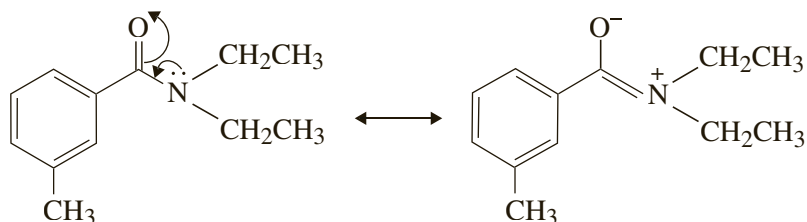

When the temperature is lowered to 0°C , the spectrum shows a pair of quartets and a pair of triplets. See the inset structures in the NMR spectrum (Figure 2) for the methylene and methyl groups, respectively.

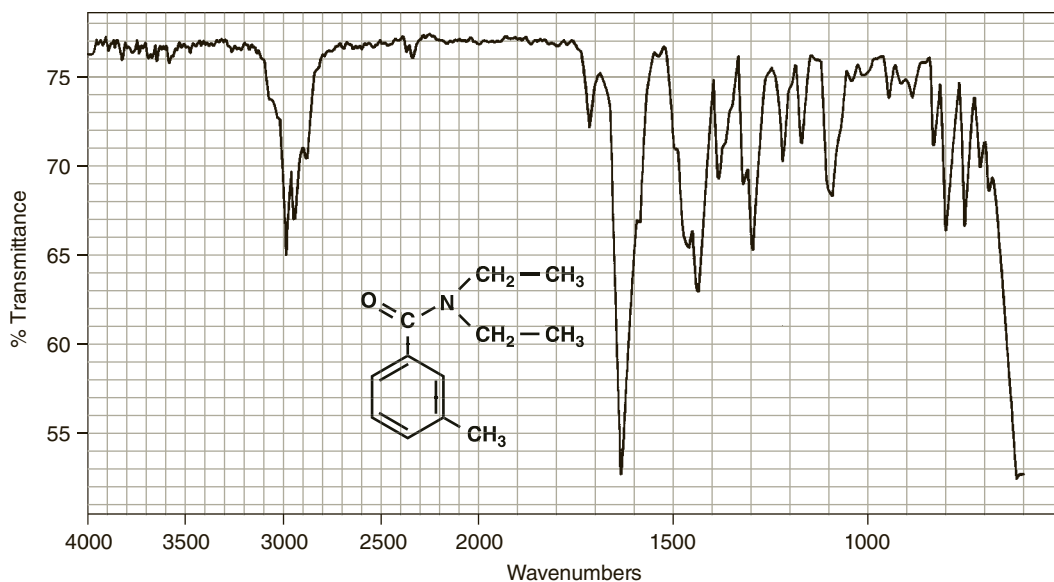


Figure 1
Infrared spectrum of *N,N*-diethyl-*m*-toluamide (neat).

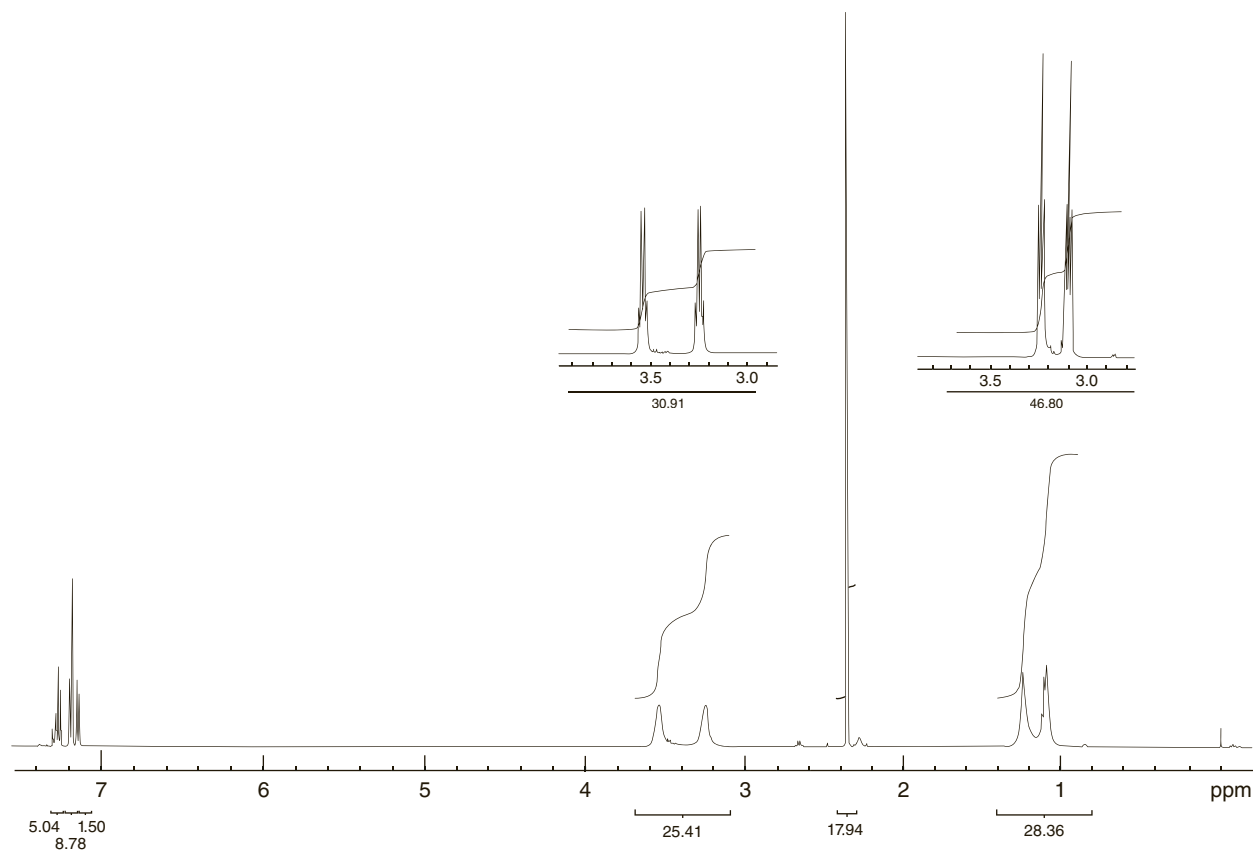


Figure 2
500 MHz ^1H NMR spectrum of *N,N*-diethyl-*m*-toluamide (CDCl_3) at 20°C (full spectrum, lower trace) and at 0°C (inset spectrum).

REFERENCE

Knoess, H. P.; Neeland, E. G. A Modified Synthesis of the Insect Repellent DEET. *J. Chem. Educ.* **1998**, *75* (Oct), 1267–1278.


QUESTIONS

1. Write an equation that describes the reaction of thionyl chloride with water.
2. What reaction would take place if the acid chloride of *m*-toluic acid were mixed with water?
3. It may be possible that some *m*-toluic acid may remain unreacted or may have formed from the hydrolysis of the acid chloride during the course of the reaction. Explain how the sodium hydroxide mixture removes unreacted carboxylic acid from the mixture. Give an equation with your answer.
4. Write a mechanism for each step in the preparation of *N,N*-diethyl-*m*-toluamide.
5. Interpret each of the principal peaks in the infrared spectrum of *N,N*-diethyl-*m*-toluamide.


Sulfa Drugs


The history of chemotherapy extends as far back as 1909 when Paul Ehrlich first used the term. Although Ehrlich's original definition of chemotherapy was limited, he is recognized as one of the giants of medicinal chemistry. **Chemotherapy** might be defined as "the treatment of disease by chemical reagents." It is preferable that these chemical reagents exhibit a toxicity toward only the pathogenic organism and not toward both the organism and the host. A chemotherapeutic agent is most useful if it does not poison the patient at the same time that it cures the patient's disease!

In 1932, the German dye manufacturing firm I. G. Farbenindustrie patented a new drug, Prontosil. Prontosil is a red azo dye, and it was first prepared for its dye properties. Remarkably, it was discovered that Prontosil showed antibacterial action when it was used to dye wool. This discovery led to studies of Prontosil as a drug capable of inhibiting the growth of bacteria. The following year, Prontosil was successfully used against staphylococcal septicemia, a blood infection. In 1935, Gerhard Domagk published the results of his research, which indicated that Prontosil was capable of curing streptococcal infections in mice and rabbits. Prontosil was shown to be active against a wide variety of bacteria in later work. This important discovery, which paved the way for a tremendous amount of research on the chemotherapy of bacterial infections, earned Domagk the 1939 Nobel Prize in medicine, but an order from Hitler prevented Domagk from accepting the honor.

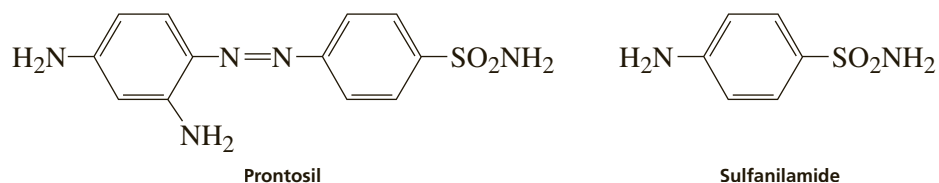


Prontosil is an effective antibacterial substance *in vivo*, that is, when injected into a living animal. Prontosil is not medicinally active when the drug is tested *in vitro*, that is, on a bacterial culture grown in the laboratory. In 1935, the research group at the Pasteur Institute in Paris headed by J. Tréfouël learned that Prontosil is metabolized in animals to sulfanilamide. Sulfanilamide had been known since 1908. Experiments with sulfanilamide showed that it had the same action as Prontosil *in vivo* and that it was also active *in vitro*, where Prontosil was known to be inactive. It was concluded that the active portion of the Prontosil molecule was the sulfanilamide moiety. This discovery led to an explosion of interest in sulfonamide derivatives. Well over a thousand sulfonamide substances were prepared within a few years of these discoveries.

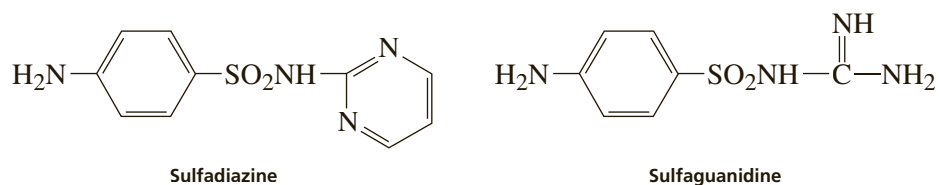

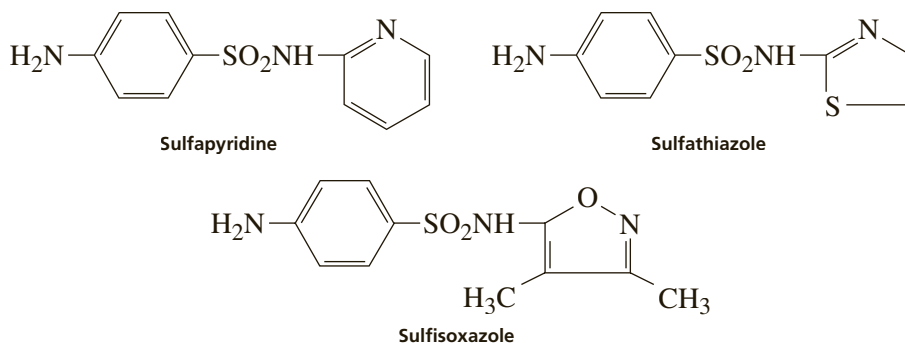

Although many sulfonamide compounds were prepared, only a relative few showed useful antibacterial properties. As the first useful antibacterial drugs, these few medicinally active sulfonamides, or **sulfa drugs**, became the wonder drugs of their day. An antibacterial drug may be either **bacteriostatic** or **bactericidal**. A bacteriostatic drug suppresses the growth of bacteria; a bactericidal drug kills bacteria. Strictly speaking, the sulfa drugs are bacteriostatic. The structures of some of the most common sulfa drugs are shown here. These more complex sulfa drugs have various important applications. Although they do not have the simple structure characteristic of sulfanilamide, they tend to be less toxic than the simpler compound.

Sulfa drugs began to lose their importance as generalized antibacterial agents when production of antibiotics in large quantity began. In 1929, Sir Alexander Fleming made his famous discovery of **penicillin**. In 1941, penicillin was first used successfully to treat humans. Since that time, the study of antibiotics has spread to molecules that bear little or no structural similarity to the sulfonamides. Besides penicillin derivatives, antibiotics that are derivatives of **tetracycline**, including Aureomycin and Terramycin, were also discovered. These newer antibiotics have high activity against bacteria, and they do not usually have the severe unpleasant side effects of many of the sulfa drugs. Nevertheless, the sulfa drugs are still widely used in treating malaria, tuberculosis, leprosy, meningitis, pneumonia, scarlet fever, plague, respiratory infections, and infections of the intestinal and urinary tracts.

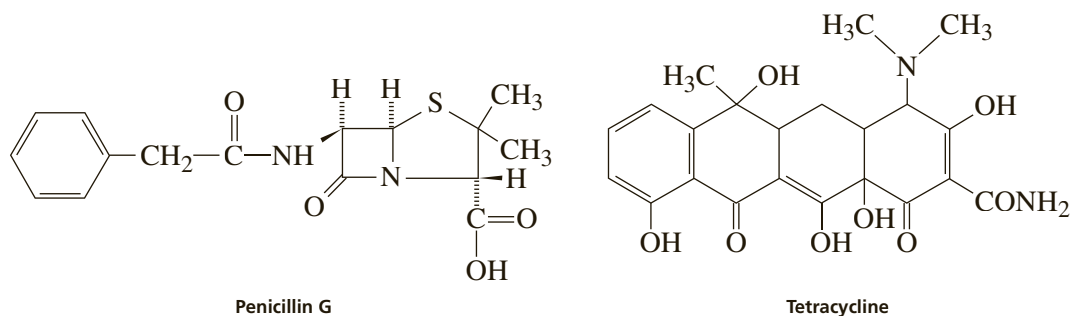


Even though the importance of sulfa drugs has declined, studies of how these materials act provide very interesting insights into how chemotherapeutic substances might behave. In 1940, Woods and Fildes discovered that *p*-aminobenzoic acid (PABA) inhibits the action of sulfanilamide. They concluded that sulfanilamide and PABA, because of their structural similarity, must compete with each other within the organism even though they cannot carry out the same chemical function. Further studies indicated that sulfanilamide does not kill bacteria, but inhibits their growth. In order to grow, bacteria require an enzyme-catalyzed reaction that

uses folic acid as a cofactor. Bacteria synthesize folic acid, using PABA as one of the components. When sulfanilamide is introduced into the bacterial cell, it competes with PABA for the active site of the enzyme that carries out the incorporation of PABA into the molecule of folic acid. Because sulfanilamide and PABA compete for an active site due to their structural similarity and because sulfanilamide cannot carry out the chemical transformations characteristic of PABA once it has formed a complex with the enzyme, sulfanilamide is called a **competitive inhibitor of the enzyme**. The enzyme, once it has formed a complex with sulfanilamide, is incapable of catalyzing the reaction required for the synthesis of folic acid. Without folic acid, the bacteria cannot synthesize the nucleic acids required for growth. As a result, bacterial growth is arrested until the body's immune system can respond and kill the bacteria.

One might well ask the question, "Why, when someone takes sulfanilamide as a drug, doesn't it inhibit the growth of *all* cells, bacterial and human alike?" The answer is simple. Animal cells cannot synthesize folic acid. Folic acid must be a part of the diet of animals and is therefore an essential vitamin. Because animal cells receive their fully synthesized folic acid molecules through the diet, only the bacterial cells are affected by the sulfanilamide and only their growth is inhibited.

For most drugs, a detailed picture of their mechanism of action is unavailable. The sulfa drugs, however, provide a rare example from which we can theorize how other therapeutic agents carry out their medicinal activity.

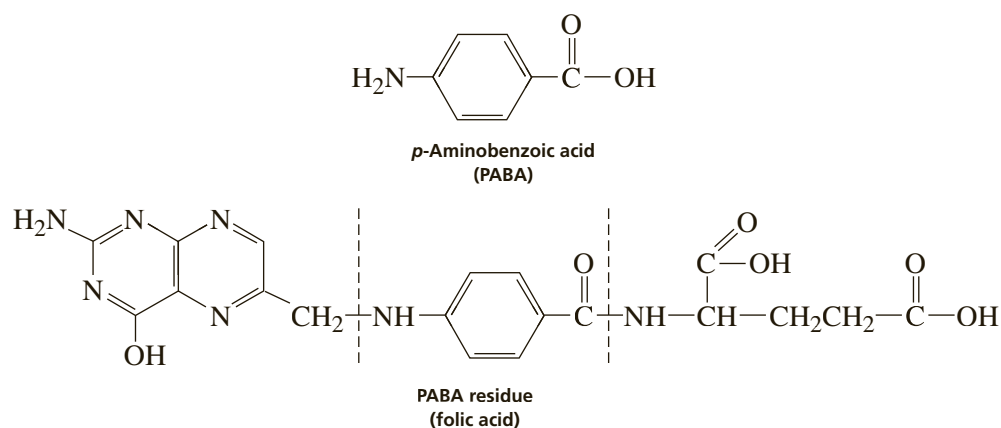

REFERENCES

- Amundsen, L. H. Sulfanilamide and Related Chemotherapeutic Agents. *J. Chem. Educ.* **1942**, *19*, 167.
- Evans, R. M. *The Chemistry of Antibiotics Used in Medicine*; Pergamon Press: London, 1965.
- Fieser, L. F.; Fieser, M. Chemotherapy. In *Topics in Organic Chemistry*; Reinhold: New York, 1963.
- Garrod, L. P.; O'Grady, F. *Antibiotics and Chemotherapy*. E. and S. Livingstone, Ltd., Edinburgh, 1968.
- Mandell, G. L.; Sande, M. A. The Sulfonamides. In *The Pharmacological Basis of Therapeutics*, 8th ed; Gilman, A. G.; Rall, T. W.; Nies, A. S.; Taylor, P., Eds.; Pergamon Press: New York, 1990.
- Sementsov, A. The Medical Heritage from Dyes. *Chemistry* **1966**, *39* (Nov), 20.
- Zahner, H.; Maas, W. K. *Biology of Antibiotics*; Springer-Verlag: Berlin, 1972.

Sulfa Drugs: Preparation of Sulfanilamide

Crystallization


Protecting groups

Testing the action of drugs on bacteria

Preparation of a sulfonamide

Aromatic substitution


In this experiment, you will prepare the sulfa drug sulfanilamide by the following synthetic scheme. The synthesis involves converting acetanilide to the intermediate *p*-acetamidobenzenesulfonyl chloride in Step 1. This intermediate is converted to sulfanilamide by way of *p*-acetamidobenzenesulfonamide in Step 2.

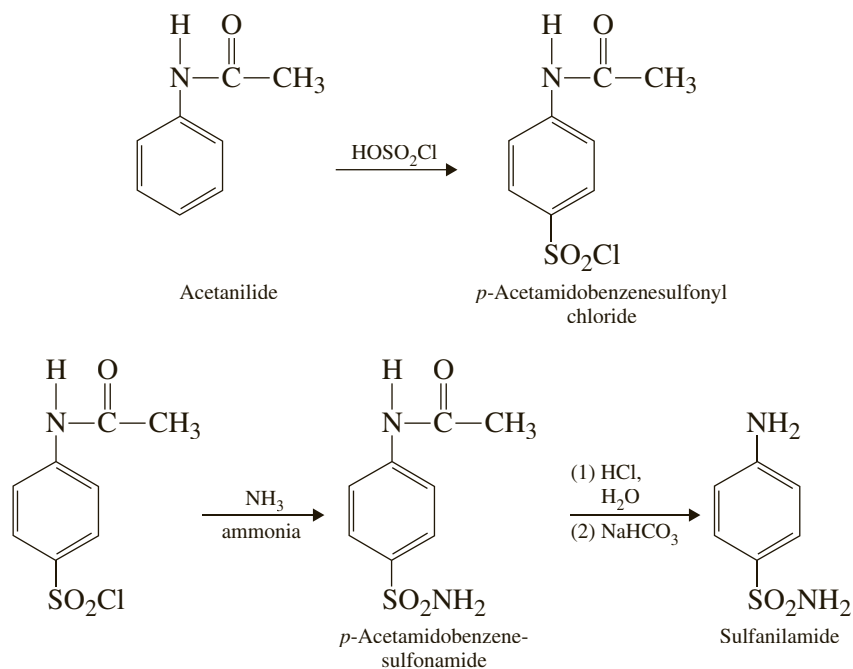

Acetanilide, which can easily be prepared from aniline, is allowed to react with chlorosulfonic acid to yield *p*-acetamidobenzenesulfonyl chloride. The acetamido group directs substitution almost totally to the *para* position. The reaction is an example of an electrophilic aromatic substitution reaction. Two problems would result if aniline itself were used in the reaction. First, the amino group in aniline would be protonated in strong acid to become a *meta* director; and, second, the chlorosulfonic acid would react with the amino group rather than with the ring, to give $\text{C}_6\text{H}_5\text{—NHSO}_3\text{H}$. For these reasons, the amino group has been “protected” by


acetylation. The acetyl group will be removed in the final step, after it is no longer needed, to regenerate the free amino group present in sulfanilamide.

p-Acetamidobenzenesulfonyl chloride is isolated by adding the reaction mixture to ice water, which decomposes the excess chlorosulfonic acid. This intermediate is fairly stable in water; nevertheless, it is converted slowly to the corresponding sulfonic acid ($\text{Ar}-\text{SO}_3\text{H}$). Thus, it should be isolated as soon as possible from the aqueous medium by filtration.

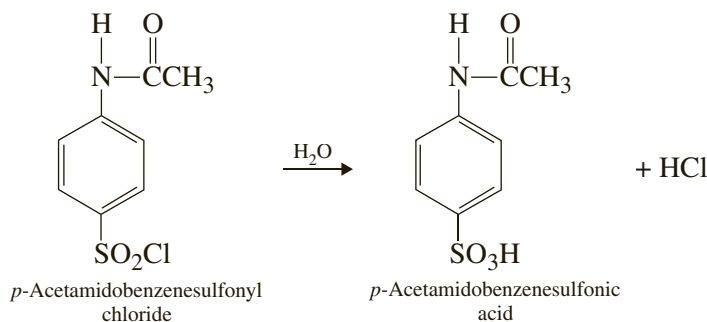

The intermediate sulfonyl chloride is converted to *p*-acetamidobenzene-sulfonamide by a reaction with aqueous ammonia (Step 2). Excess ammonia neutralizes the hydrogen chloride produced. The only side reaction is the hydrolysis of the sulfonyl chloride to *p*-acetamidobenzenesulfonic acid.

The protecting acetyl group is removed by acid-catalyzed hydrolysis to generate the hydrochloride salt of the product, sulfanilamide. Note that of the two amide linkages present, only the carboxylic acid amide (acetamido group) was cleaved, not the sulfonic acid amide (sulfonamide). The salt of the sulfa drug is converted to sulfanilamide when the base, sodium bicarbonate, is added.

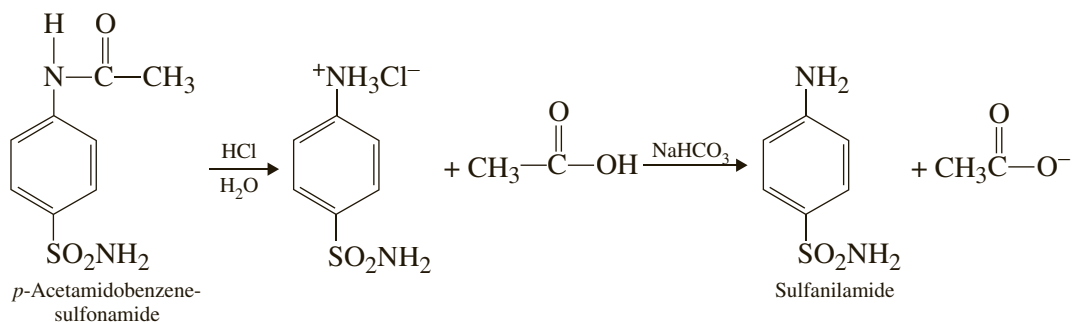

REQUIRED READING

Review:	Technique 7	Reaction Methods
	Technique 8	Filtration, Sections 8.3 and 8.7
	Technique 11	Crystallization, Section 11.4
	Technique 25	Infrared Spectroscopy, Sections 25.4 and 25.5
New:	Essay	Sulfa Drugs

SPECIAL INSTRUCTIONS

Chlorosulfonic acid must be handled with care because it is a corrosive liquid and reacts violently with water. The *p*-acetamidobenzenesulfonyl chloride should be used during the same laboratory period in which it is prepared. It is unstable and will not survive long storage. The sulfa drug may be tested on several kinds of bacteria (Instructor's Manual).

SUGGESTED WASTE DISPOSAL

Aqueous filtrates should be placed in the container provided for this purpose. Place organic wastes in the nonhalogenated waste container.

PROCEDURE

Part A. *p*-Acetamidobenzenesulfonyl Chloride

The Reaction Apparatus

Assemble the apparatus as shown in Technique 7, Figure 7.6A (inset) using dry glassware. You will need a 5-mL conical vial, an air condenser, and a drying tube, which will be used as a gas trap. Prepare the drying tube for use as a gas trap by packing the tube loosely with dry glass wool (Technique 7, Section 7.8A). Moisten the glass wool slightly with several drops of water. The moistened glass wool traps the hydrogen chloride that is evolved in the reaction. Attach the 5-mL conical vial after the acetanilide and chlorosulfonic acid have been added, as directed in the following paragraph. You should adjust the temperature of the aluminum block to about 110°C for use later in the experiment.

Reaction of Acetanilide with Chlorosulfonic Acid

Place 0.18 g of acetanilide in the dry 5-mL conical vial and connect the air condenser but not the drying tube. Melt the acetanilide (mp 113°C) by heating the vial in a community sand bath or aluminum block set to about 160°C. Remove the vial from the heating source and swirl the heavy oil while holding the vial at an angle so that it is deposited uniformly on the cone-shaped bottom of the vial. Allow the conical vial to cool to room temperature and then cool it further in an ice-water bath. (Don't place the hot vial directly into the ice-water bath without prior cooling, or the vial will crack.)

CAUTION



Chlorosulfonic acid is an extremely noxious and corrosive chemical and should be handled with care. Use only dry glassware with this reagent. Should the chlorosulfonic acid be spilled on your skin, wash it off immediately with water. Be very careful when washing any glassware that has come in contact with chlorosulfonic acid. Even a small amount of the acid will react vigorously with water and may splatter. Wear safety glasses.

Remove the air condenser. In a hood, transfer 0.50 mL of chlorosulfonic acid ClSO_2OH ($MW = 116.5$, $d = 1.77 \text{ g/mL}$) to the acetanilide in the conical vial using the graduated pipette provided. Reattach the air condenser and drying tube. Allow the mixture to stand for 5 minutes and then heat the reaction vial in the aluminum

block at about 110°C for 10 minutes to complete the reaction. Remove the vial from the aluminum block. Allow the vial to cool to the touch and then cool it in an ice-water bath.

Isolation of p-Acetamidobenzenesulfonyl Chloride

The operations described in this paragraph should be conducted as rapidly as possible because the *p*-acetamidobenzenesulfonyl chloride reacts with water. Add 3 g of crushed ice to a 20-mL beaker. In a hood, transfer the cooled reaction mixture dropwise (it may splatter somewhat) with a Pasteur pipette onto the ice while stirring the mixture with a glass stirring rod. (The remaining operations in this paragraph may be completed at your laboratory bench.) Rinse the conical vial with a few drops of cold water and transfer the contents to the beaker containing the ice. Stir the precipitate to break up the lumps and then filter the *p*-acetamidobenzenesulfonyl chloride on a Hirsch funnel (Technique 8, Section 8.3, and Figure 8.5). Rinse the conical vial and beaker with two 1-mL portions of ice water. Use the rinse water to wash the crude product on the funnel. Any remaining solid in the conical vial should be left there because this vial is used again in the next section. Do not stop here. Convert the solid into *p*-acetamidobenzenesulfonamide in the same laboratory period.

Part B. Sulfanilamide *Preparation of p-Acetamidobenzenesulfonamide*

Prepare a hot water bath at 70°C. Place the crude *p*-acetamidobenzenesulfonyl chloride into the original 5-mL conical vial and add 1.1 mL of dilute ammonium hydroxide solution.¹ Stir the mixture well with a spatula and reattach the air condenser and drying tube (gas trap) using fresh, moistened glass wool. Heat the mixture in the hot water bath for 10 minutes. Allow the conical vial to cool to the touch and place it in an ice-water bath for several minutes. Collect the *p*-acetamidobenzenesulfonamide on a Hirsch funnel and rinse the vial and product with a small amount of ice water. You may stop here.

Hydrolysis of p-Acetamidobenzenesulfonamide


Transfer the solid into the conical vial and add 0.53 mL of dilute hydrochloric acid solution.² Attach the air condenser and heat the mixture in an aluminum block at about 130°C until all the solid has dissolved. Then heat the solution for an additional 5 minutes. Allow the mixture to cool to room temperature. If a solid (unreacted starting material) appears, heat the mixture for several minutes at 130°C. When the vial has cooled to room temperature, no further solids should appear.

Isolation of Sulfanilamide

With a Pasteur pipette, transfer the solution to a 20-mL beaker. While stirring with a glass rod, cautiously add dropwise a slurry of 0.5 g of sodium bicarbonate in about 1 mL of water to the mixture in the beaker. Foaming will occur after each addition of the bicarbonate solution because of carbon dioxide evolution. Allow gas evolution to cease before making the next addition. Eventually, sulfanilamide will begin to precipitate. At this point, begin to check the pH of the solution. Add the aqueous sodium bicarbonate until the pH of the solution is between 4 and 6. Cool the mixture thoroughly in an ice-water bath. Collect the sulfanilamide on a Hirsch funnel and rinse the beaker and solid with about 0.5 mL of cold water. Allow the solid to air dry on the Hirsch funnel for several minutes using suction.

¹ Prepared by mixing 11.0 mL of concentrated ammonium hydroxide with 11.0 mL of water.

² Prepared by mixing 7.0 mL of water with 3.6 mL of concentrated hydrochloric acid.

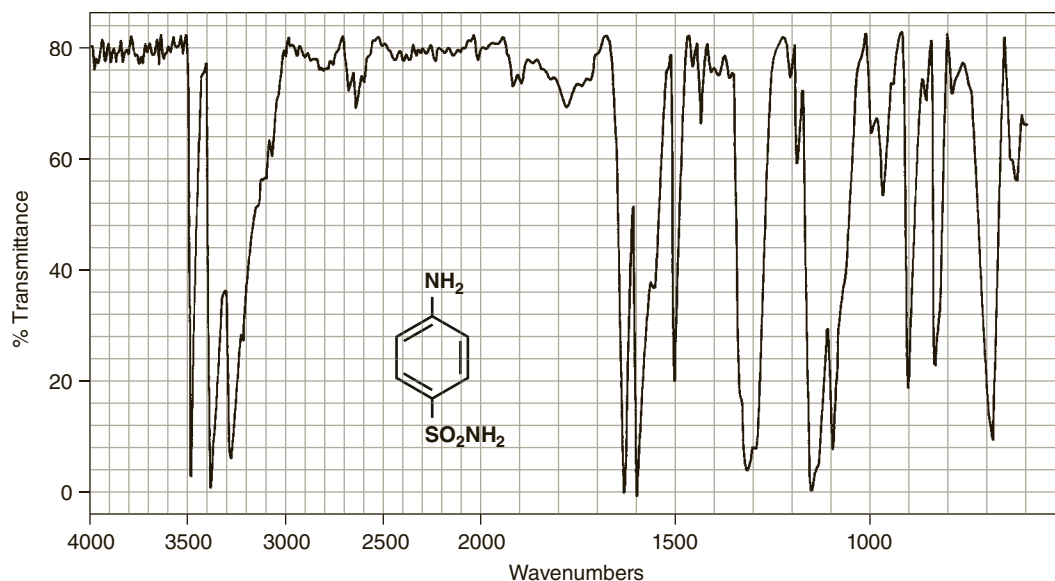

Infrared spectrum of sulfanilamide, KBr.

Crystallization of Sulfanilamide


Weigh the crude product and crystallize it from hot water (use 1.0 to 1.2 mL water/0.1 g) using a Craig tube (Technique 11, Section 11.4, and Figure 11.6). Step 2 in Figure 11.6 (removal of insoluble impurities) should not be required in this crystallization. Let the purified product dry until the next laboratory period.

Yield Calculation, Melting Point, and Infrared Spectrum

Weigh the dry sulfanilamide and calculate the percentage yield ($MW = 172.2$). Determine the melting point (pure sulfanilamide melts at $163\text{--}164^\circ\text{C}$). At the option of the instructor, obtain the infrared spectrum using the dry film method (Technique 25, Section 25.4) or by an alternative method (Technique 25, Section 25.5). Compare your infrared spectrum with the one reproduced here. Submit the sulfanilamide to the instructor in a labeled vial or save it for the tests with bacteria (see Instructor's Manual).

QUESTIONS

1. Write an equation showing how excess chlorosulfonic acid is decomposed in water.
2. In the preparation of sulfanilamide, why was aqueous sodium bicarbonate, rather than aqueous sodium hydroxide, used to neutralize the solution in the final step?
3. At first glance, it might seem possible to prepare sulfanilamide from sulfanilic acid by the set of reactions shown here.

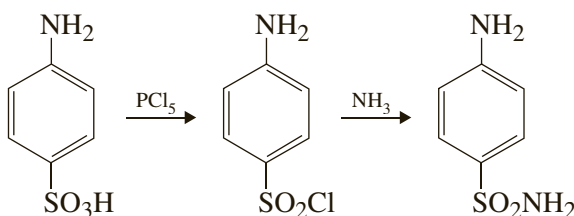

When the reaction is conducted in this way, however, a polymeric product is produced after Step 1. What is the structure of the polymer? Why does *p*-acetamidobenzoyl chloride not produce a polymer?


Polymers and Plastics

Chemically, plastics are composed of chainlike molecules of high molecular weight called **polymers**. Polymers have been built up from simpler chemicals called **monomers**. The word *poly* is defined as “many,” *mono* means “one,” and *mer* indicates “units.” Thus, many monomers are combined to give a polymer. A different monomer or combination of monomers is used to manufacture each type or family of polymers. There are two broad classes of polymers: addition and condensation. Both types are described here.


Many polymers (plastics) produced in the past were of such low quality that they gained a bad reputation. The plastics industry now produces high-quality materials that are increasingly replacing metals in many applications. They are used in many products such as clothes, toys, furniture, machine components, paints, boats, automobile parts, and even artificial organs. In the automobile industry, metals have been replaced with plastics to help reduce the overall weight of cars and to help reduce corrosion. This reduction in weight helps improve gas mileage. Epoxy resins can even replace metal in engine parts.

CHEMICAL STRUCTURES OF POLYMERS

Basically, a polymer is made up of many repeating molecular units formed by sequential addition of monomer molecules to one another. Many monomer molecules of A, say 1,000 to 1 million, can be linked to form a gigantic polymeric molecule:

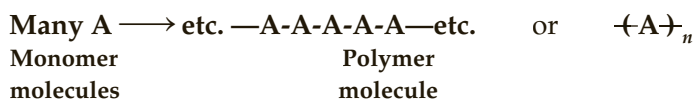

Monomers that are different can also be linked to form a polymer with an alternating structure. This type of polymer is called a copolymer.

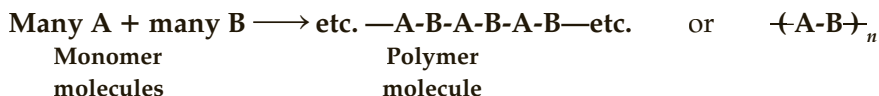


TYPES OF POLYMERS

For convenience, chemists classify polymers in several main groups, depending on the method of synthesis.

1. **Addition polymers** are formed by a reaction in which monomer units simply add to one another to form a long-chain (generally linear or branched) polymer. The monomers usually contain carbon-carbon double bonds. Examples of synthetic addition polymers include polystyrene (Styrofoam),


polytetrafluoroethylene (Teflon), polyethylene, polypropylene, polyacrylonitrile (Orlon, Acrilan, Creslan), poly(vinyl chloride) (PVC), and poly(methyl methacrylate) (Lucite, Plexiglas). The process can be represented as follows:

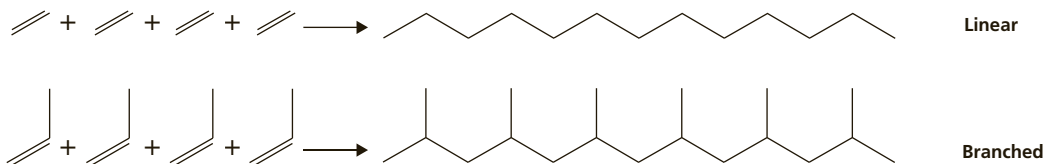

2. **Condensation polymers** are formed by the reaction of bifunctional or polyfunctional molecules, with the elimination of some small molecule (such as water, ammonia, or hydrogen chloride) as a by-product. Familiar examples of synthetic condensation polymers include polyesters (Dacron, Mylar), polyamides (nylon), polyurethanes, and epoxy resin. Natural condensation polymers include polyamino acids (protein), cellulose, and starch. The process can be represented as follows:



3. **Cross-linked polymers** are formed when long chains are linked in one gigantic, three-dimensional structure with tremendous rigidity. Addition and condensation polymers can exist with a cross-linked network, depending on the monomers used in the synthesis. Familiar examples of cross-linked polymers are Bakelite, rubber, and casting (boat) resin. The process can be represented as follows:

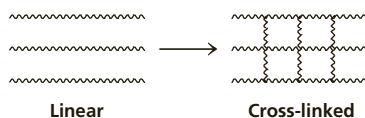

Linear and crossed-linked polymers.


THERMAL CLASSIFICATION OF POLYMERS

Industrialists and technologists often classify polymers as either thermoplastics or thermoset plastics rather than as addition or condensation polymers. This classification takes into account their thermal properties.

1. **Thermal properties of thermoplastics.** Most addition polymers and many condensation polymers can be softened (melted) by heat and reformed (molded) into other shapes. Industrialists and technologists often refer to these types of polymers as **thermoplastics**. Weaker, noncovalent bonds (dipole-dipole and London dispersion) are broken during the heating. Technically, thermoplastics are the materials we call plastics. Thermoplastics may be repeatedly melted and recast into new shapes. They may be recycled as long as degradation does not occur during reprocessing.

Some addition polymers, such as poly(vinyl chloride), are difficult to melt and process. Liquids with high boiling points, such as dibutyl phthalate, are added to the polymer to separate the chains from each other. These compounds are called **plasticizers**. In effect, they act as lubricants that neutralize the attractions that exist between chains. As a result, the polymer can be melted at a lower temperature to aid in processing. In addition, the polymer becomes more flexible

at room temperature. By varying the amount of plasticizer, poly(vinyl chloride) can range from a very flexible, rubberlike material to a very hard substance.

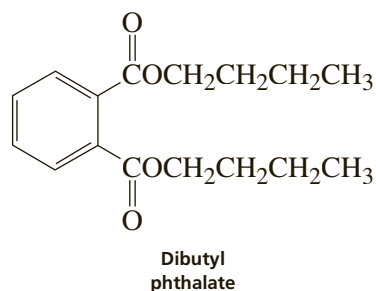


Phthalate plasticizers are volatile compounds of low molecular weight. Part of the new car smell comes from the odor of these materials as they evaporate from the vinyl upholstery. The vapor often condenses on the windshield as an oily film. After some time, the vinyl material may lose enough plasticizer to cause it to crack.

2. Thermal properties of thermoset plastics. Industrialists use the term **thermoset** plastics to describe materials that melt initially but on further heating become permanently hardened. Once formed, thermoset materials cannot be softened and remolded without destruction of the polymer, because covalent bonds are broken. Thermoset plastics cannot be recycled. Chemically, thermoset plastics are cross-linked polymers. They are formed when long chains are linked in one gigantic, three-dimensional structure with tremendous rigidity.

Polymers can also be classified in other ways; for example, many varieties of rubber are often referred to as *elastomers*, Dacron is a *fiber*, and poly(vinyl acetate) is an *adhesive*. The addition and condensation classifications are used in this essay.

ADDITION POLYMERS

By volume, most of the polymers prepared in industry are of the addition type. The monomers generally contain a carbon-carbon double bond. The most important example of an addition polymer is the well-known polyethylene, for which the monomer is ethylene. Countless numbers (n) of ethylene molecules are linked in long-chain polymeric molecules by breaking the pi bond and creating two new single bonds between the monomer units. The number of recurring units may be large or small, depending on the polymerization conditions.

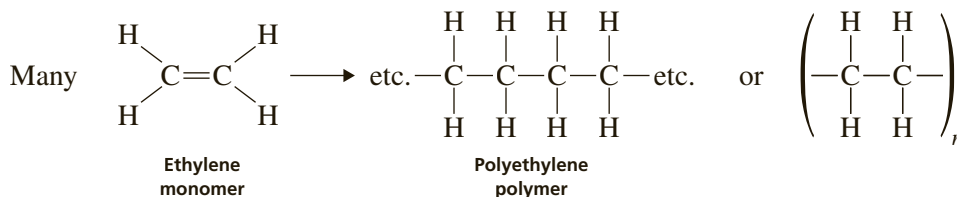


This reaction can be promoted by heat, pressure, and a chemical catalyst. The molecules produced in a typical reaction vary in the number of carbon atoms in their chains. In other words, a mixture of polymers of varying length, rather than a pure compound, is produced.

Polyethylenes with linear structures can pack together easily and are referred to as high-density polyethylenes. They are fairly rigid materials. Low-density

polyethylenes consist of branched-chain molecules, with some cross-linking in the chains. They are more flexible than the high-density polyethylenes. The reaction conditions and the catalysts that produce polyethylenes of low and high density are quite different. The monomer, however, is the same in each case.

Another example of an addition polymer is polypropylene. In this case, the monomer is propylene. The polymer that results has a branched methyl on alternate carbon atoms of the chain.

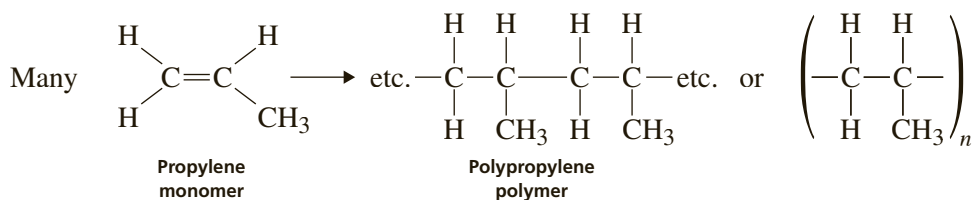


A number of common addition polymers are shown in Table 1. Some of their principal uses are also listed. The last three entries in the table all have a carbon-carbon double bond remaining after the polymer is formed. These bonds activate or participate in a further reaction to form cross-linked polymers called *elastomers*; this term is almost synonymous with *rubber*, because elastomers are materials with common characteristics.

CONDENSATION POLYMERS

Condensation polymers, for which the monomers contain more than one type of functional group, are more complex than addition polymers. In addition, most condensation polymers are copolymers made from more than one type of monomer. Recall that addition polymers, in contrast, are all prepared from substituted ethylene molecules. The single functional group in each case is one or more double bonds, and a single type of monomer is generally used.

Dacron, a polyester, can be prepared by causing a dicarboxylic acid to react with a bifunctional alcohol (a diol):

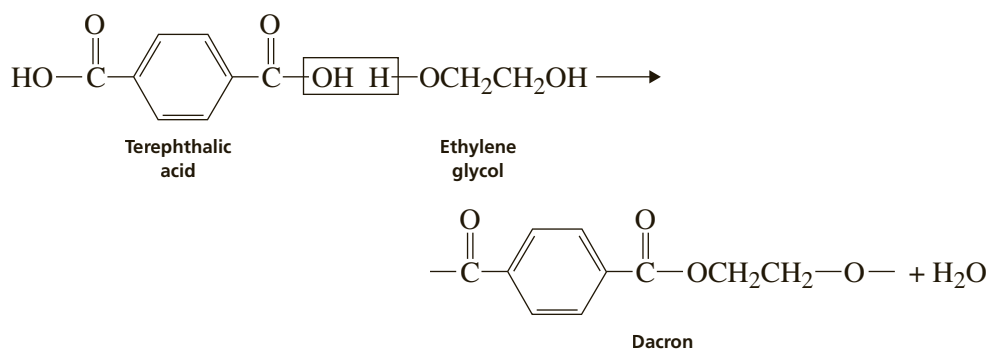

Nylon 6-6, a polyamide, can be prepared by causing a dicarboxylic acid to react with a bifunctional amine.

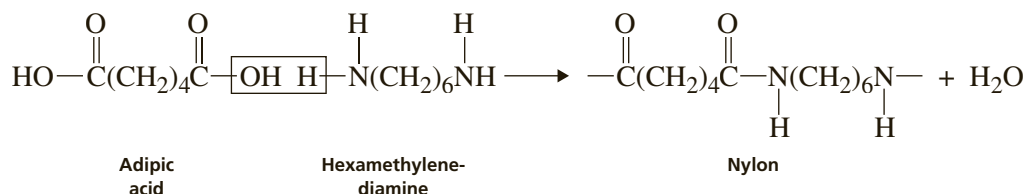


Table 1 Addition Polymers

Example	Monomer(s)	Polymer	Uses
Polyethylene	$\text{CH}_2=\text{CH}_2$	$-\text{CH}_2-\text{CH}_2-$	Most common and important polymer; bags, insulation for wires, squeeze bottles
Polypropylene	$\begin{array}{c} \text{CH}_2 = \text{CH} \\ \\ \text{CH}_3 \end{array}$	$-\text{CH}_2-\underset{\text{CH}_3}{\text{CH}}-$	Fibers, indoor-outdoor carpets, bottles
Polystyrene	$\begin{array}{c} \text{CH}_2 = \text{CH} \\ \\ \text{C}_6\text{H}_5 \end{array}$	$-\text{CH}_2-\underset{\text{C}_6\text{H}_5}{\text{CH}}-$	Styrofoam, inexpensive household goods, inexpensive molded objects
Poly(vinyl chloride) (PVC)	$\begin{array}{c} \text{CH}_2 = \text{CH} \\ \\ \text{Cl} \end{array}$	$-\text{CH}_2-\underset{\text{Cl}}{\text{CH}}-$	Synthetic leather, clear bottles, floor covering, phonograph records, water pipe
Polytetrafluoroethylene (Teflon)	$\text{CF}_2=\text{CF}_2$	$-\text{CF}_2-\text{CF}_2-$	Nonstick surfaces, chemically resistant films
Poly(methyl methacrylate) (Lucite, Plexiglas)	$\begin{array}{c} \text{CO}_2\text{CH}_3 \\ \\ \text{CH}_2 = \text{C} \\ \\ \text{CH}_3 \end{array}$	$-\text{CH}_2-\underset{\text{CH}_3}{\overset{\text{CO}_2\text{CH}_3}{\text{C}}}-$	Unbreakable "glass," latex paints
Polyacrylonitrile (Orlon, Acrilan, Creslan)	$\begin{array}{c} \text{CH}_2 = \text{CH} \\ \\ \text{CN} \end{array}$	$-\text{CH}_2-\underset{\text{CN}}{\text{CH}}-$	Fiber used in sweaters, blankets, carpets
Poly(vinyl acetate) (PVA)	$\begin{array}{c} \text{CH}_2 = \text{CH} \\ \\ \text{OCCH}_3 \\ \\ \text{O} \end{array}$	$-\text{CH}_2-\underset{\text{OCCH}_3}{\text{CH}}-$	Adhesives, latex paints, chewing gum, textile coatings
Natural rubber	$\begin{array}{c} \text{CH}_3 \\ \\ \text{CH}_2 = \text{CCH} = \text{CH}_2 \end{array}$	$-\text{CH}_2-\underset{\text{CH}_3}{\text{C}} = \text{CH}-\text{CH}_2-$	Polymer cross-linked with sulfur (vulcanization)
Polychloroprene (neoprene rubber)	$\begin{array}{c} \text{Cl} \\ \\ \text{CH}_2 = \text{CCH} = \text{CH}_2 \end{array}$	$-\text{CH}_2-\underset{\text{Cl}}{\text{C}} = \text{CH}-\text{CH}_2-$	Cross-linked with ZnO; resistant to oil and gasoline
Styrene butadiene rubber (SBR)	$\begin{array}{c} \text{CH}_2 = \text{CH} \\ \\ \text{C}_6\text{H}_5 \\ \text{CH}_2 = \text{CHCH} = \text{CH}_2 \end{array}$	$-\text{CH}_2\underset{\text{C}_6\text{H}_5}{\text{CH}}-\text{CH}_2\text{CH} = \text{CHCH}_2-$	Cross-linked with peroxides; most common rubber, used for tires; 25% styrene, 75% butadiene


Notice, in each case, that a small molecule, water, is eliminated as a product of the reaction. Several other condensation polymers are listed in Table 2. Linear (or branched) chain polymers, as well as cross-linked polymers, are produced in condensation reactions.

Table 2 Condensation Polymers

Example	Monomers	Polymer	Uses
Polyamides (nylon)	$\begin{array}{c} \text{O} \quad \text{O} \\ \parallel \quad \parallel \\ \text{HO}(\text{C}(\text{H}_2)_n\text{C}(\text{O})\text{OH} \\ \text{H}_2\text{N}(\text{C}(\text{H}_2)_n\text{NH}_2 \end{array}$	$-\text{C}(\text{H}_2)_n\text{C}(\text{O})-\text{NH}(\text{C}(\text{H}_2)_n\text{NH}-$	Fibers, molded objects
Polyesters (Dacron, Mylar, Fortrel)	$\begin{array}{c} \text{O} \quad \text{O} \\ \parallel \quad \parallel \\ \text{HO}(\text{C}_6\text{H}_4\text{C}(\text{O})\text{COH} \\ \text{HO}(\text{C}(\text{H}_2)_n\text{OH} \end{array}$	$-\text{C}(\text{H}_2)_n\text{C}(\text{O})-\text{O}(\text{C}_6\text{H}_4\text{C}(\text{O})\text{CO}-$	Linear polyesters, fibers, recording tape
Polyesters (Glyptal resin)	$\begin{array}{c} \text{O} \quad \text{O} \\ \parallel \quad \parallel \\ \text{HO}(\text{C}_6\text{H}_3\text{C}(\text{O})\text{OC}(\text{O})\text{C}_6\text{H}_3\text{C}(\text{O})\text{OH} \\ \text{HOCH}_2\text{CH}(\text{OH})\text{CH}_2\text{OH} \end{array}$	$-\text{C}_6\text{H}_3\text{C}(\text{O})\text{OC}(\text{O})\text{C}_6\text{H}_3\text{C}(\text{O})\text{OCH}_2\text{CH}(\text{OH})\text{CH}_2\text{O}-$	Cross-linked polyester, paints
Polyesters (casting resin)	$\begin{array}{c} \text{O} \quad \text{O} \\ \parallel \quad \parallel \\ \text{HO}(\text{C}(\text{H}_2)_n\text{C}(\text{O})\text{CH}=\text{CH}(\text{C}(\text{O})\text{OH} \\ \text{HO}(\text{C}(\text{H}_2)_n\text{OH} \end{array}$	$-\text{C}(\text{H}_2)_n\text{C}(\text{O})\text{CH}=\text{CH}(\text{C}(\text{O})\text{O}(\text{C}(\text{H}_2)_n\text{O}-$	Cross-linked with styrene and peroxide; fiberglass boat resin
Phenol-formaldehyde resin (Bakelite)	$\begin{array}{c} \text{OH} \\ \\ \text{C}_6\text{H}_5 \\ \\ \text{CH}_2=\text{O} \end{array}$	$-\text{CH}_2-\text{C}_6\text{H}_3(\text{OH})_2-\text{CH}_2-\text{C}_6\text{H}_3(\text{OH})_2-\text{CH}_2-$	Mixed with fillers; molded electrical goods, adhesives, laminates, varnishes
Cellulose acetate*	$\begin{array}{c} \text{CH}_2\text{OH} \\ \\ \text{C}_6\text{H}_4 \\ \\ \text{OH} \\ \\ \text{OH} \\ \\ \text{CH}_3\text{COOH} \end{array}$	$\begin{array}{c} \text{CH}_2\text{OAc} \\ \\ \text{C}_6\text{H}_4 \\ \\ \text{OAc} \\ \\ \text{OAc} \end{array}$	Photographic film
Silicones	$\begin{array}{c} \text{CH}_3 \\ \\ \text{Cl}-\text{Si}-\text{Cl} \quad \text{H}_2\text{O} \\ \\ \text{CH}_3 \end{array}$	$-\text{O}-\text{Si}(\text{CH}_3)_2-\text{O}-$	Water-repellent coatings, temperature-resistant fluids and rubbers (CH_3SiCl_3 cross-links in water)
Polyurethanes	$\begin{array}{c} \text{CH}_3 \\ \\ \text{C}_6\text{H}_4 \\ \\ \text{N}=\text{C}=\text{O} \\ \\ \text{N}=\text{C}=\text{O} \\ \\ \text{HO}(\text{C}(\text{H}_2)_n\text{OH} \end{array}$	$\begin{array}{c} \text{CH}_3 \\ \\ \text{NHC}(\text{O})-\text{O}(\text{C}(\text{H}_2)_n\text{O}- \\ \\ \text{NHC}(\text{O})-\text{O}(\text{C}(\text{H}_2)_n\text{O}- \\ \\ \text{O} \end{array}$	Rigid and flexible foams, fibers

*Cellulose, a polymer of glucose, is used as the monomer.

The nylon structure contains the amide linkage at regular intervals:

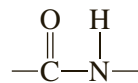


This type of linkage is extremely important in nature because of its presence in proteins and polypeptides. Proteins are gigantic polymeric substances made up of monomer units of amino acids. They are linked by the peptide (amide) bond.

Other important natural condensation polymers are starch and cellulose. They are polymeric materials made up of the sugar monomer glucose. Another important natural condensation polymer is the DNA molecule. A DNA molecule is made up of the sugar deoxyribose linked with phosphates to form the backbone of the molecule.

Polycarbonates are another important type of condensation polymer widely used in the marketplace. Since they are a thermoplastic material, they can be easily molded into a number of different products. Polycarbonates have outstanding high-impact resistance, which make them ideal for use as “unbreakable” water bottles and food storage containers. They also have outstanding optical properties, which make them highly desirable for lenses in high-impact eyewear. Since polycarbonates have low scratch resistance, a hard coating is usually applied to the surface of lenses. Polycarbonates have replaced glass in many applications because of their durability, clarity, breakage resistance, and light weight. Polycarbonates share some characteristics with the older, more-established material, poly(methyl methacrylate); the structure of this material is shown in Table 1. However, polycarbonates are stronger and more durable than poly(methyl methacrylate). Although more expensive, polycarbonates can be identified by looking for the number 7 stamped on the bottoms of containers. Category 7 is the catch-all code for “other” plastics (see Table 3).



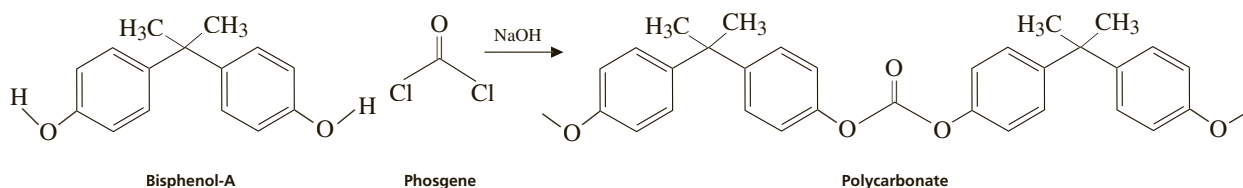




The most common type of polycarbonate is made from bisphenol-A (BPA). One way of preparing this plastic involves the reaction of bisphenol-A and phosgene in the presence of sodium hydroxide.




Bisphenol-A is very much in the news today. There is fear that some of this monomer may find its way into food. The major concern is possible contamination from baby bottles made of polycarbonate. The worry is that bisphenol-A may be formed from the break-down of polycarbonate used in baby bottles. If this happens, then bisphenol-A would contaminate infant formula or milk in the bottles and be ingested by babies. In the laboratory setting, bisphenol-A also appears to be released from animal cages made from waste polycarbonate. It appears when water leaches small amounts of it out of the plastic. The study suggests that bisphenol-A may be responsible for enlargement of the reproductive organs of female mice. In the past, these studies have been disputed by the chemical industry, which argued that the average dose of bisphenol-A is far too low to be harmful—a finding initially supported by the Federal Drug Administration (FDA).

Recent animal studies have suggested, however, that even small doses of bisphenol-A exposure can cause a number of health risks and may mimic the female hormone estrogen. The study suggests that feminizing effects can develop in

Table 3 Code System for Plastic Materials

Code	Polymer	Uses
 PETE	Poly(ethylene terephthalate) (PET) $-O-CH_2-CH_2-O-C(=O)-C_6H_4-C(=O)-$	Soft-drink bottles
 HDPE	High-density polyethylene $-CH_2-CH_2-CH_2-CH_2-$	Milk and beverage containers, products in squeeze bottle
 V	Vinyl/poly(vinyl chloride) (PVC) $-CH_2-\underset{\text{Cl}}{\text{CH}}-CH_2-\underset{\text{Cl}}{\text{CH}}-$	Some shampoo containers, bottles with cleaning materials in them
 LDPE	Low-density polyethylene $-CH_2-CH_2-CH_2-CH_2-$ with some branches	Thin plastic bags, some plastic wrap
 PP	Polypropylene $-CH_2-\underset{\text{CH}_3}{\text{CH}}-CH_2-\underset{\text{CH}_3}{\text{CH}}-$	Heavy-duty, microwaveable containers used in kitchens
 PS	Polystyrene $-CH_2-\underset{\text{C}_6\text{H}_5}{\text{CH}}-CH_2-\underset{\text{C}_6\text{H}_5}{\text{CH}}-$	Beverage/foam cups, window in envelopes
 Other	All other resins, layered multimaterials, containers made of different materials	Some ketchup bottles, snack packs, mixtures where top differs from bottom

fetuses and infants. Studies reported in the *Journal of the American Medical Society* found that higher levels of bisphenol-A in adults were associated with greater incidences of diabetes and cardiovascular problems. In October 2008, the FDA found its original assessment to be flawed. In the meantime most manufacturers of water bottles have changed their formulation. On April 18, 2008, Health Canada announced that bisphenol-A is “toxic to human health.” Canada is the first country to make this designation. Eastman’s Triton® was accepted as a suitable alternative in August 2008 by Health Canada. This material is described as a “copolyester” by the manufacturer. The alcohol components in the Triton polyesters are often mixtures of 2,2,4,4-tetramethylcyclobutane-1,4-diol and 1,4-cyclohexanedimethanol. Often the dicarboxylic acid component is terephthalic acid. Other manufacturers may use some 1,3-propanediol in this polyester formulations, along with tetramethylcyclobutanediol.

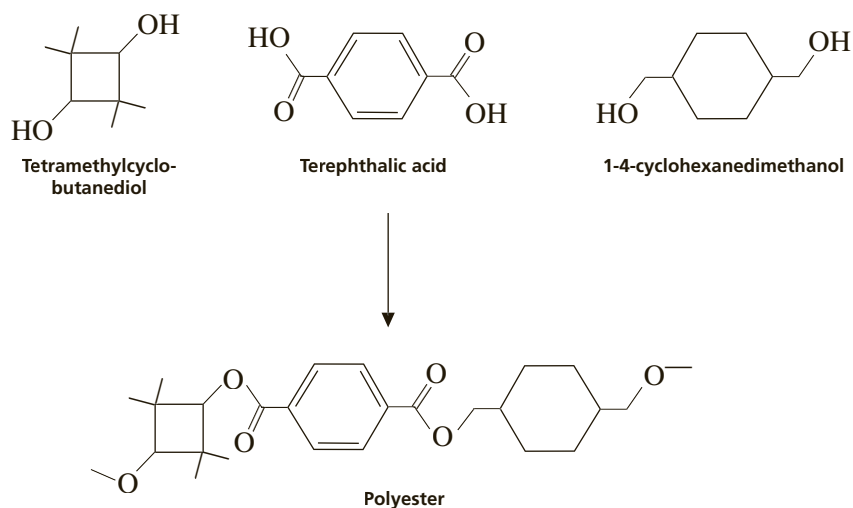

Unfortunately, bisphenol-A (BPA) is also one of the components in the most common type of epoxy resin. BPA-based epoxy resins are often applied to the inside of food and soft-drink cans in order to form a protective coating between the metal can and the food material inside. It turns out that the epoxy resins adhere readily to metal containers, and other potential substitutes do not appear to be as suitable for food-contact purposes. The Food and Drug Administration has not recommended discontinuing the use of BPA-based epoxy, at least at this time.

Ideally, we should either recycle all our wastes or not produce the waste in the first place. Plastic waste consists of about 55% polyethylene and polypropylene, 20% polystyrene, and 11% PVC. All these polymers are thermoplastics and can be recycled. They can be resoftened and remolded into new goods. Unfortunately, thermosetting plastics (crosslinked polymers) cannot be remelted. They decompose on high-temperature heating. Thus, thermosetting plastics should not be used for “disposable” purposes. To recycle plastics effectively, we must sort the materials according to the various types. The plastics industry has introduced a code system consisting of seven categories for the common plastics used in packaging. The code is conveniently stamped on the bottom of the container. Using these codes, consumers can separate the plastics into groups for recycling purposes. These codes are listed in Table 3, together with the most common uses around the home. Notice that the seventh category is a miscellaneous one, called “Other.”

It is quite amazing that so few different plastics are used in packaging. The most common ones are polyethylene (low and high density), polypropylene, polystyrene, and poly(ethylene terephthalate). All of these materials can easily be recycled because they are thermoplastics. Incidentally, vinyls (polyvinyl chloride) are becoming less common in packaging.

REFERENCES

- Ainsworth, S. J. Plastics Additives. *Chem. Eng. News* **1992**, (Aug 31), 34–55.
- Burfield, D. R. Polymer Glass Transition Temperatures. *J. Chem. Educ.* **1987**, *64*, 875.
- Carraher, C. E., Jr.; Hess, G.; Sperling, L. H. Polymer Nomenclature—or What’s in a Name? *J. Chem. Educ.* **1987**, *64*, 36.
- Carraher, C. E., Jr.; Seymour, R. B. Physical Aspects of Polymer Structure: A Dictionary of Terms. *J. Chem. Educ.* **1986**, *63*, 418.
- Carraher, C. E., Jr.; Seymour, R. B. Polymer Properties and Testing—Definitions. *J. Chem. Educ.* **1987**, *64*, 866.

- Carraher, C. E., Jr.; Seymour, R. B. Polymer Structure—Organic Aspects (Definitions). *J. Chem. Educ.* **1988**, *65*, 314.
- Erickson, B. E.; Ritter, S. K.; Tullo, A. H. Pressure on Plasticizers. *Chem. Eng. News* **2015**, *93* (15), (June 22), 10–20.
- Fried, J. R. The Polymers of Commercial Plastics. *Plast. Eng.* **1982**, (Jun), 49–55.
- Fried, J. R. Polymer Properties in the Solid State. *Plast. Eng.* **1982**, (Jul), 27–37.
- Fried, J. R. Molecular Weight and Its Relation to Properties. *Plast. Eng.* **1982**, (Aug), 27–33.
- Fried, J. R. Elastomers and Thermosets. *Plast. Eng.* **1983**, (Mar), 67–73.
- Fried, J. R.; Yeh, E. B. Polymers and Computer Alchemy. *Chemtech* **1993**, *23*, (Mar), 35–40.
- Goodall, B. L. The History and Current State of the Art of Propylene Polymerization Catalysts. *J. Chem. Educ.* **1986**, *63*, 191.
- Harris, F. W.; et al. State of the Art: Polymer Chemistry. *J. Chem. Educ.* **1981**, *58*, (Nov). (This issue contains 17 papers on polymer chemistry. The series covers structures, properties, mechanisms of formation, methods of preparation, stereochemistry, molecular weight distribution, rheological behavior of polymer melts, mechanical properties, rubber elasticity, block and graft copolymers, organometallic polymers, fibers, ionic polymers, and polymer compatibility.)
- Hileman, B. Bisphenol A May Trigger Human Breast Cancer: Study in Rats Provides Strongest Case Yet Against Common Environmental Chemical. *Chem. Eng. News* [Online] **2006**, *84* (Dec 6), <http://pubs.acs.org/cen/news/84/i50/8450bisphenol.html>
- Howdeshell, K. L.; Peterman, P. H.; Judy, B. M.; Taylor, J. A.; Orazio, C. E.; Ruhlen, R. L.; Vom Saal, F. S.; Welshons, W. V. Bisphenol A Is Released from Used Polycarbonate Animal Cages into Water at Room Temperature. *Environ. Health Perspect.* **2003**, (Jul), 1180–1187.
- Jacoby, M. Trading Places with Bisphenol A, BPA-based Polymers' Useful Properties Make Them Tough Materials to Replace. *Chem. Eng. News* **2008**, *86* (50), (Dec 15), 31–33.
- Jordan, R. F. Cationic Metal–Alkyl Olefin Polymerization Catalysts. *J. Chem. Educ.* **1988**, *65*, 285.
- Kauffman, G. B. Wallace Hume Carothers and Nylon, the First Completely Synthetic Fiber. *J. Chem. Educ.* **1988**, *65*, 803.
- Kauffman, G. B. Rayon: The First Semi-Synthetic Fiber Product. *J. Chem. Educ.* **1993**, *70*, 887.
- Kauffman, G. B.; Seymour, R. B. Elastomers I. Natural Rubber. *J. Chem. Educ.* **1990**, *67*, 422.
- Kauffman, G. B.; Seymour, R. B. Elastomers II. Synthetic Rubbers. *J. Chem. Educ.* **1991**, *68*, 217.
- Lang, I. A.; Galloway, T. S.; Scarlett, A.; Henley, W. E.; Depledge, M.; Wallace, R. B.; Melzer, D. Association of Urinary Bisphenol A Concentration with Medical Disorders and Laboratory Abnormalities in Adults. *J. Am. Med. Assoc.* **300** (1), 1303–1310.
- Morse, P. M. New Catalysts Renew Polyolefins. *Chem. Eng. News* **1998**, *76*, (Jul 6), 11–16.
- Seymour, R. B. Polymers Are Everywhere. *J. Chem. Educ.* **1988**, *65*, 327.
- Seymour, R. B. Alkenes and Their Derivatives: The Alchemists' Dream Come True. *J. Chem. Educ.* **1989**, *66*, 670.
- Seymour, R. B.; Kauffman, G. B. Polymer Blends: Superior Products from Inferior Materials. *J. Chem. Educ.* **1992**, *69*, 646.
- Seymour, R. B.; Kauffman, G. B. Polyurethanes: A Class of Modern Versatile Materials. *J. Chem. Educ.* **1992**, *69*, 909.
- Seymour, R. B.; Kauffman, G. B. The Rise and Fall of Celluloid. *J. Chem. Educ.* **1992**, *69*, 311.
- Seymour, R. B.; Kauffman, G. B. Thermoplastic Elastomers. *J. Chem. Educ.* **1992**, *69*, 967.
- Stevens, M. P. Polymer Additives: Chemical and Aesthetic Property Modifiers. *J. Chem. Educ.* **1993**, *70*, 535.
- Stevens, M. P. Polymer Additives: Mechanical Property Modifiers. *J. Chem. Educ.* **1993**, *70*, 444.
- Stevens, M. P. Polymer Additives: Surface Property and Processing Modifiers. *J. Chem. Educ.* **1993**, *70*, 713.
- Thayer, A. M. Metallocene Catalysts Initiate New Era in Polymer Synthesis. *Chem. Eng. News* **1995**, *73* (Sept 11), 15–20.
- Waller, F. J. Fluoropolymers. *J. Chem. Educ.* **1989**, *66*, 487.
- Webster, O. W. Living Polymerization Methods. *Science* **1991**, *251*, 887.

Preparation and Properties of Polymers: Polyester, Nylon, and Polystyrene

Condensation polymers

Addition polymers

Cross-linked polymers

Infrared spectroscopy

In this experiment, the syntheses of two polyesters (Experiment 52A), nylon (Experiment 52B), and polystyrene (Experiment 52C) will be described. These polymers represent important commercial plastics. They also represent the main classes of polymers: condensation (linear polyester, nylon), addition (polystyrene), and cross-linked (Glyptal polyester). Infrared spectroscopy is used in Experiment 52D to determine the structure of polymers.

REQUIRED READING

Review: Technique 25 Infrared Spectroscopy, Section 25B

New: Essay Polymers and Plastics

SPECIAL INSTRUCTIONS

Experiments 52A, 52B, and 52C all involve toxic vapors. Each experiment should be conducted in a well-ventilated hood. The styrene used in Experiment 52C irritates the skin and eyes. Avoid breathing its vapors. Styrene must be dispensed and stored in a hood. Benzoyl peroxide is flammable and may detonate on impact or on heating.


SUGGESTED WASTE DISPOSAL

The test tubes containing the polyester polymers from Experiment 52A should be placed in a box designated for disposal of these samples. The nylon from Experiment 52B should be washed thoroughly with water and placed in a waste container. The liquid wastes from Experiment 52B (nylon) should be poured into a container designated for disposal of these wastes. The polystyrene prepared in Experiment 52C should be placed in the container designated for solid wastes.

52A EXPERIMENT 52A


Polyesters


Linear and cross-linked polyesters will be prepared in this experiment. Both are examples of condensation polymers. The linear polyester is prepared as follows:

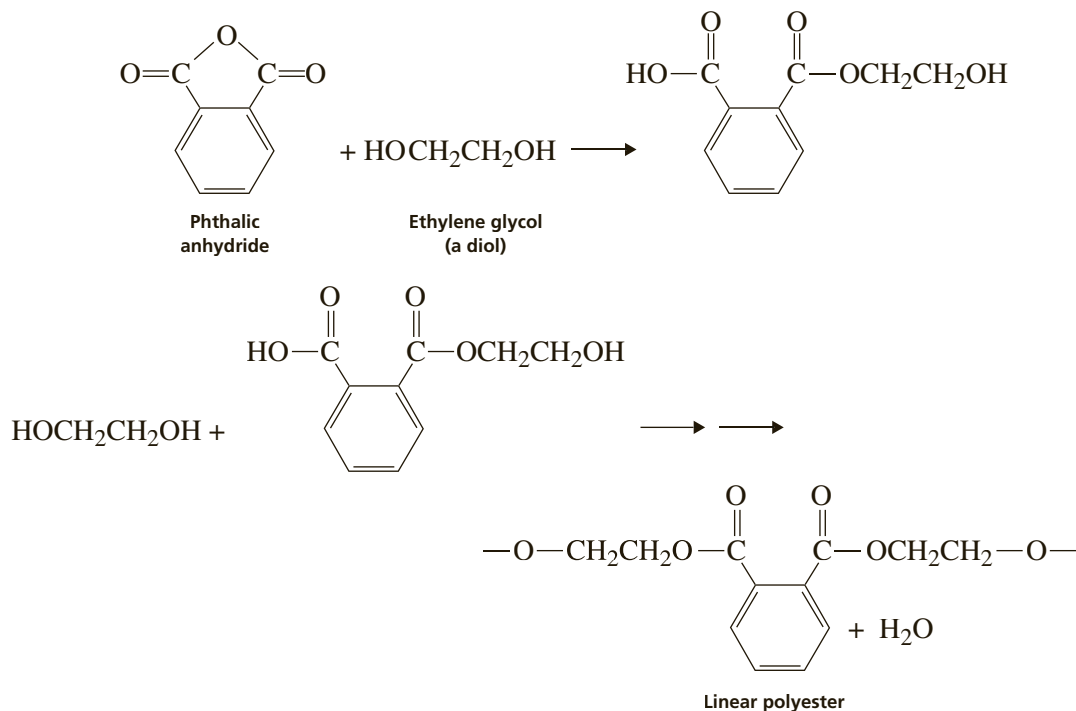

This linear polyester is isomeric with Dacron, which is prepared from terephthalic acid and ethylene glycol (see the preceding essay). Dacron and the linear polyester made in this experiment are both thermoplastics.

If more than two functional groups are present in one of the monomers, the polymer chains can be linked to one another (cross-linked) to form a three-dimensional network. Such structures are usually more rigid than linear structures and are useful in making paints and coatings. They may be classified as thermosetting plastics. The polyester Glyptal is prepared as follows:

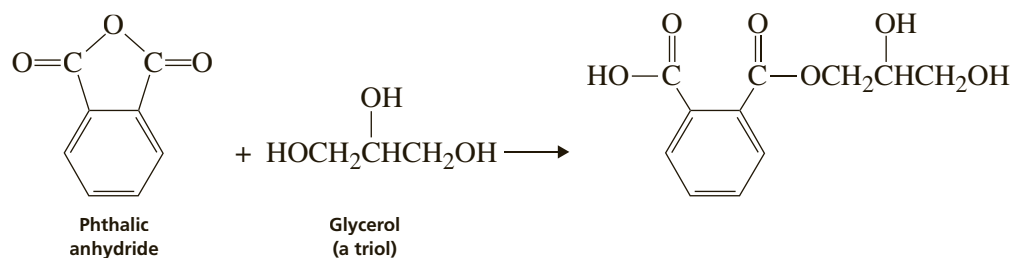

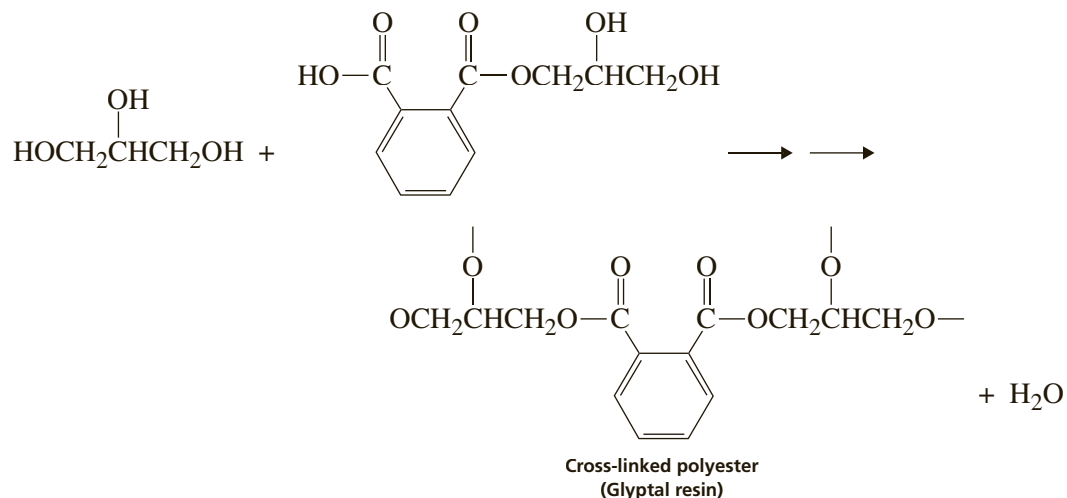

The reaction of phthalic anhydride with a diol (ethylene glycol) is described in the procedure. This linear polyester is compared with the cross-linked polyester (Glyptal) prepared from phthalic anhydride and a triol (glycerol).

PROCEDURE

Place 1 g of phthalic anhydride and 0.05 g of sodium acetate in each of two test tubes. To one tube, add 0.4 mL of ethylene glycol and to the other, add 0.4 mL of glycerol. Clamp both tubes so that they can be heated simultaneously with a flame. Heat the tubes gently until the solutions appear to boil (water is eliminated during the esterification); then continue heating for 5 minutes.

If you are performing the optional infrared analysis of the polymer, immediately save a sample of the polymer formed from ethylene glycol only. After removing a sample for infrared spectroscopy, allow the two test tubes to cool and compare the viscosity and brittleness of the two polymers. The test tubes cannot be cleaned.

Optional Exercise: Infrared Spectroscopy



Lightly coat a watch glass with stopcock grease. Pour some of the *hot* polymer from the tube containing ethylene glycol; use a wooden applicator stick to spread the polymer on the surface to create a thin film of the polymer. Remove the polymer from the watch glass and save it for Experiment 52D.

52B EXPERIMENT 52B

Polyamide (Nylon)

Reaction of a dicarboxylic acid, or one of its derivatives, with a diamine leads to a linear polyamide through a condensation reaction. Commercially, nylon 6-6 (so called because each monomer has six carbons) is made from adipic acid and hexamethylenediamine. In this experiment, you will use the acid chloride instead of adipic acid. The acid chloride is dissolved in cyclohexane, and this is added *carefully* to hexamethylenediamine dissolved in water. These liquids do not mix, so two layers will form. The polymer can then be drawn out continuously to form a long

strand of nylon. Imagine how many molecules have been linked in this long strand! It is a fantastic number.

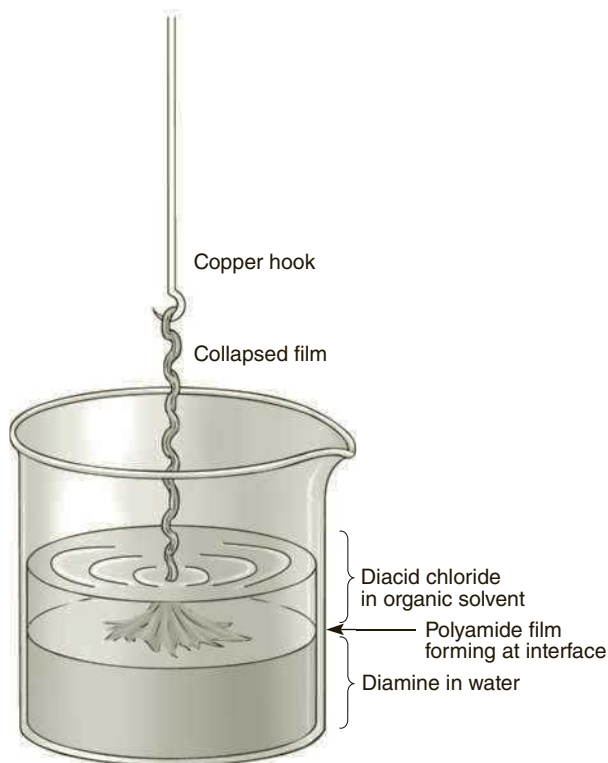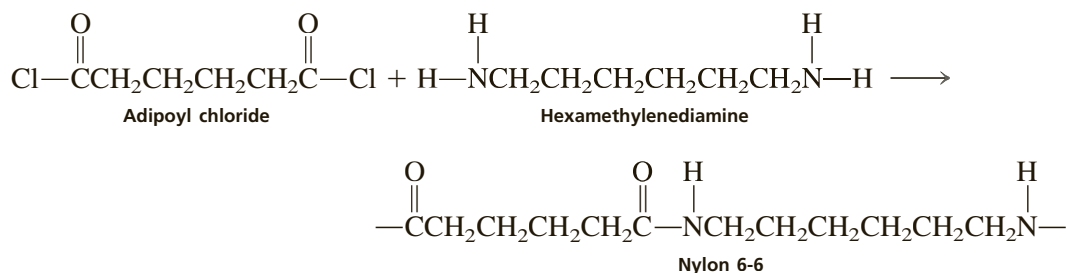

Preparation of nylon.

PROCEDURE


Pour 10 mL of a 5% aqueous solution of hexamethylenediamine (1,6-hexanediamine) into a 50-mL beaker. Add 10 drops of 20% sodium hydroxide solution. Carefully add 10 mL of a 5% solution of adipoyl chloride in cyclohexane to the solution by pouring it down the wall of the slightly tilted beaker. Two layers will form (see figure), and there will be an immediate formation of a polymer film at the liquid–liquid interface. Using a copper-wire hook (a 6-inch piece of wire bent at one end), gently free the walls of the beaker from polymer strings. Then hook the mass at the center and slowly raise the wire so that polyamide forms continuously, producing a rope that can be drawn out for many feet. The strand can be broken by pulling it faster. Rinse the rope several times with water and lay it on a paper towel to dry. With the piece of wire, vigorously stir the remainder of the two-phase system to form additional polymer. Decant the liquid and wash the polymer thoroughly with water. Allow the polymer to dry. Do not discard the nylon in the sink; use a waste container.

52C EXPERIMENT 52C


Polystyrene

An addition polymer, polystyrene, will be prepared in this experiment. Reaction can be brought about by free-radical, cationic, or anionic catalysts (initiators), the first of these being the most common. In this experiment, polystyrene is prepared by free-radical-initiated polymerization.


The reaction is initiated by a free-radical source. The initiator will be benzoyl peroxide, a relatively unstable molecule, which at 80–90°C decomposes with homolytic cleavage of the oxygen–oxygen bond:

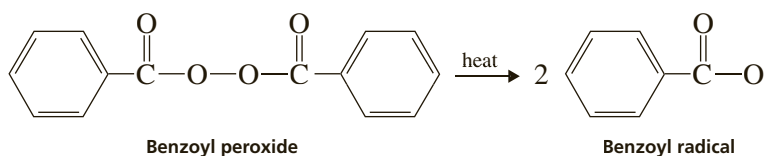

If an unsaturated monomer is present, the radical adds to it, initiating a chain reaction by producing a new free radical. If we let R stand for the initiator radical, the reaction with styrene can be represented as

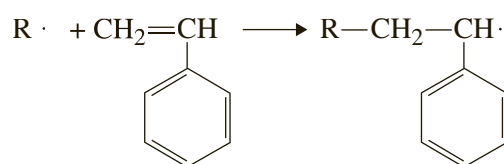

The chain continues to grow:

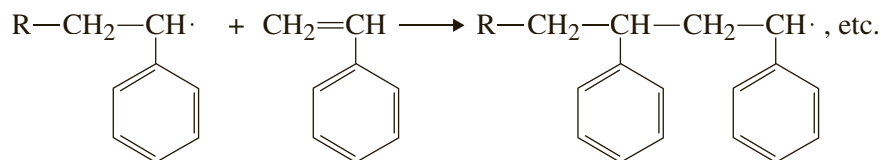

The chain can be terminated by causing two radicals to combine (either both polymer radicals or one polymer radical and one initiator radical) or by causing a hydrogen atom to become abstracted from another molecule.

PROCEDURE

Because it is difficult to clean the glassware, this experiment is best performed by the laboratory instructor. One large batch of polystyrene should be made for the entire class (at least 10 times the amounts given). After the polystyrene is prepared, a small amount will be dispensed to each student. The students will provide their own watch glass for this purpose. Perform the experiment in a hood. Place several thicknesses of newspaper in the hood.

CAUTION



Styrene vapor is very irritating to eyes, mucous membranes, and upper respiratory tract. Do not breathe the vapor and do not get it on your skin. Exposure can cause nausea and headaches. All operations with styrene must be conducted in a hood.

Benzoyl peroxide is flammable and may detonate on impact or on heating (or grinding). It should be weighed on glassine (glazed, not ordinary) paper. Clean all spills with water. Wash the glassine paper with water before discarding it.

Place 12–15 mL of styrene monomer in a 100-mL beaker and add 0.35 g of benzoyl peroxide. Heat the mixture on a hot plate until the mixture turns yellow. When the color disappears and bubbles begin to appear, immediately take the beaker of styrene off the hot plate because the reaction is exothermic (use tongs or an insulated glove). After the reaction subsides, put the beaker of styrene back on the hot plate and continue heating it until the liquid becomes very syrupy. With a stirring rod, draw out a long filament of material from the beaker. If this filament can be cleanly snapped after a few seconds of cooling, the polystyrene is ready to be poured. If the filament does not break, continue heating the mixture and repeat this process until the filament breaks easily.

If you are performing the optional infrared analysis of the polymer, immediately save a sample of the polymer. After removing a sample for infrared spectroscopy, pour the remainder of the syrupy liquid on a watch glass that has been lightly coated with stopcock grease. After being cooled, the polystyrene can be lifted from the glass surface by gently prying with a spatula.

Optional Exercise: Infrared Spectroscopy

Pour a small amount of the *hot* polymer from the beaker onto a warm watch glass (no grease) and spread the polymer with a wooden applicator stick to create a thin film of the polymer. Peel the polymer from the watch glass and save it for 52D.

52D EXPERIMENT 52D

Infrared Spectra of Polymer Samples

Infrared spectroscopy is an excellent technique for determining the structure of a polymer. For example, polyethylene and polypropylene have relatively simple spectra because they are saturated hydrocarbons. Polyesters have stretching frequencies associated with the C=O and C—O groups in the polymer chain. Polyamides (nylon) show absorptions that are characteristic for the C=O stretch and N—H stretch. Polystyrene has characteristic features of a monosubstituted aromatic compound (see Technique 25, Figure 25.12). You may determine the infrared spectra of the linear polyester from Experiment 52A and polystyrene from Experiment 52C in this part of the experiment. Your instructor may ask you to analyze a sample that you bring to the laboratory or one supplied to you.

PROCEDURE

Mounting the Samples

Prepare cardboard mounts for your polymer samples. Cut 3 × 5-inch index cards so that they fit into the sample cell holder of your infrared spectrometer. Then cut a 0.5-inch-wide × 1-inch-high rectangular hole in the center of the cardstock. Attach a polymer sample on the cardboard mount with tape.

Choices of Polymer Samples

If you have completed Experiments 52A and 52C, you can obtain the spectra of your polyester or polystyrene. Alternatively, your instructor may provide you with known or unknown polymer samples for you to analyze.

Your instructor may ask you to bring a polymer sample of your own choice. If possible, these samples should be clear and as thin as possible (similar to the thickness of plastic sandwich wrap). Good choices of plastic materials include windows from envelopes, plastic sandwich wrap, sandwich bags, soft-drink bottles, milk containers, shampoo bottles, candy wrappers, and shrink-wrap. If necessary, the samples can be heated in an oven and stretched to obtain thinner samples. If you are bringing a sample cut from a plastic container, obtain the recycling code from the bottom of the container, if one is given.

Running the Infrared Spectrum

Insert the cardboard mount into the cell holder in the spectrometer so that your polymer sample is centered in the infrared beam of the instrument. Find the thinnest place in your polymer sample. Determine the infrared spectrum of your sample. Because of the thickness of your polymer sample, many absorptions are so strong that you will not be able to see individual bands. To obtain a better spectrum, try moving the sample to a new position in the beam and rerun the spectrum. If your instrument has an ATR accessory available, you may find that you will obtain better results using this device (Technique 25, Section 25.5).

Analyzing the Infrared Spectrum


You can use the essay "Polymers and Plastics" and Technique 25 with your spectrum to help determine the structure of the polymer. Most likely, the polymers will consist of plastic materials listed in Table Three of the essay. This table lists the recycling codes for a number of household plastics used in packaging. Submit the infrared spectrum along with the structure of the polymer to your instructor. Do your spectrum and structure agree with the recycling code? Label the spectrum with the important absorption bands consistent with the structure of the polymer.

Using a Polymer Library


If your particular instrument has a polymer library, you can search the library for a match. Do this after you have made a preliminary "educated guess" as to the structure of the polymer. The library search should help confirm the structure you determined.

QUESTIONS


1. Ethylene dichloride ($\text{ClCH}_2\text{CH}_2\text{Cl}$) and sodium polysulfide (Na_2S_4) react to form a chemically resistant rubber, Thiokol A. Draw the structure of the rubber.
2. Draw the structure for the polymer produced from the monomer vinylidene chloride ($\text{CH}_2=\text{CCl}_2$).
3. Draw the structure of the copolymer produced from vinyl acetate and vinyl chloride. This copolymer is employed in some paints, adhesives, and paper coatings.

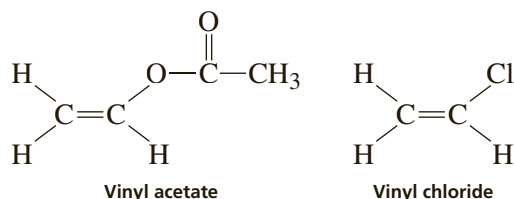

4. Isobutylene, $\text{CH}_2=\text{C}(\text{CH}_3)_2$, is used to prepare cold-flow rubber. Draw a structure for the addition polymer formed from this alkene.
5. Kel-F is an addition polymer with the following partial structure. What is the monomer used to prepare it?

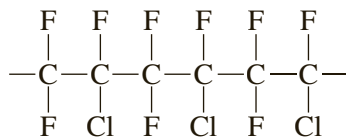


6. Maleic anhydride reacts with ethylene glycol to produce an alkyd resin. Draw the structure of the condensation polymer produced.

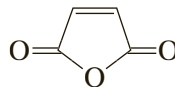


Maleic anhydride

7. Kodel is a condensation polymer made from terephthalic acid and 1,4-cyclohexanedimethanol. Write the structure of the resulting polymer.

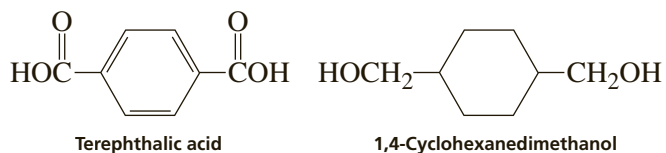

Diels–Alder Reaction and Insecticides


Since the 1930s, it has been known that the addition of an unsaturated molecule across a diene system forms a substituted cyclohexene. The original research dealing with this type of reaction was performed by Otto Diels and Kurt Alder in Germany, and the reaction is known today as the **Diels–Alder reaction**. The Diels–Alder reaction is the reaction of a **diene** with a species capable of reacting with the diene, the **dienophile**.

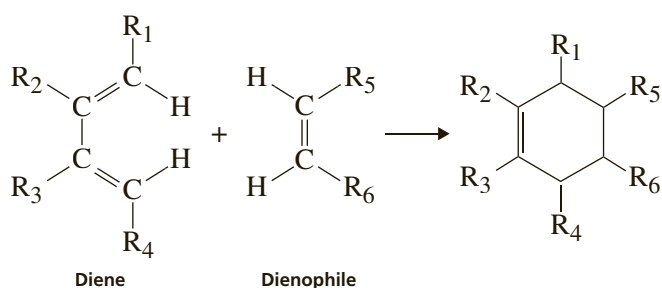


The product of the Diels–Alder reaction is usually a structure that contains a cyclohexene ring system. If the substituents as shown are simply alkyl groups or hydrogen atoms, the reaction proceeds only under extreme conditions of temperature and pressure. With more complex substituents, however, the Diels–Alder reaction may proceed at low temperatures and under mild conditions. The reaction of cyclopentadiene with maleic anhydride (Experiment 53) is an example of a Diels–Alder reaction carried out under reasonably mild conditions.

In the past, a commercially important use of the Diels–Alder reaction involved the use of hexachlorocyclopentadiene as the diene. Depending on the dienophile, a variety of chlorine-containing addition products may be synthesized. Nearly all these products were powerful insecticides. Three insecticides synthesized by the Diels–Alder reaction are shown below.

Dieldrin and Aldrin are named after Diels and Alder. These insecticides were once used against the insect pests of fruits, vegetables, and cotton; against soil insects, termites, and moths; and in the treatment of seeds. Chlordane was used in veterinary medicine against insect pests of animals, including fleas, ticks, and lice. These insecticides are seldom used today.

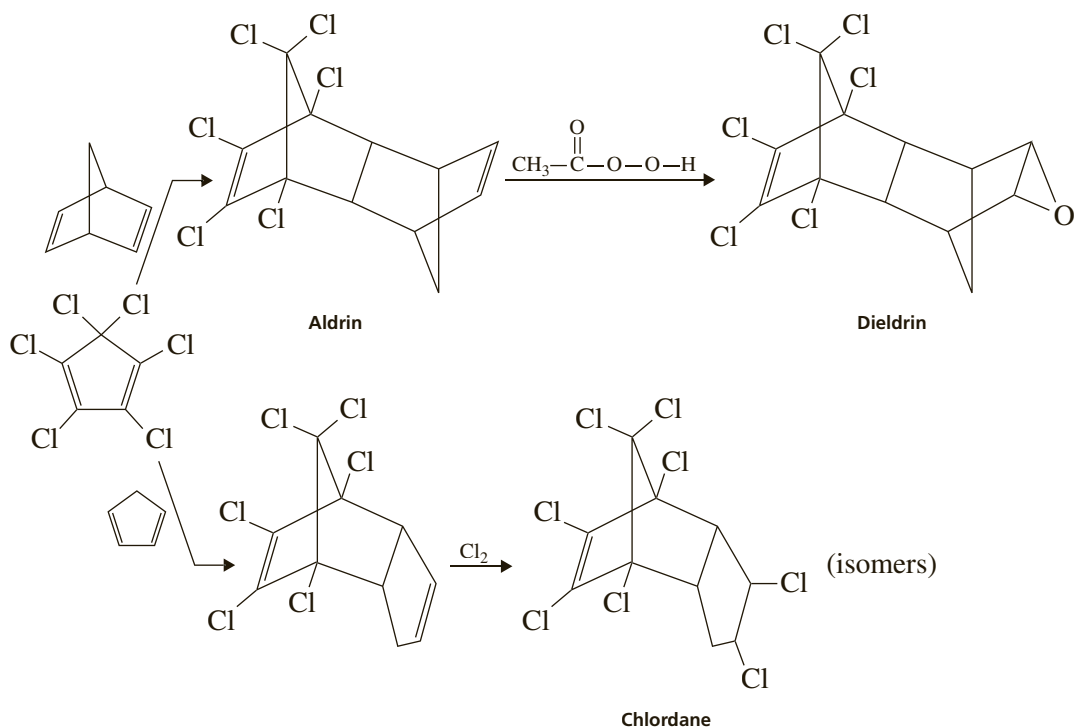

The best known insecticide, DDT, is not prepared by the Diels–Alder reaction, but is nevertheless the best illustration of the difficulties that were experienced when chlorinated insecticides were used indiscriminately. DDT was first synthesized in 1874, and its insecticidal properties were first demonstrated in 1939. It is easily synthesized commercially, with inexpensive reagents.

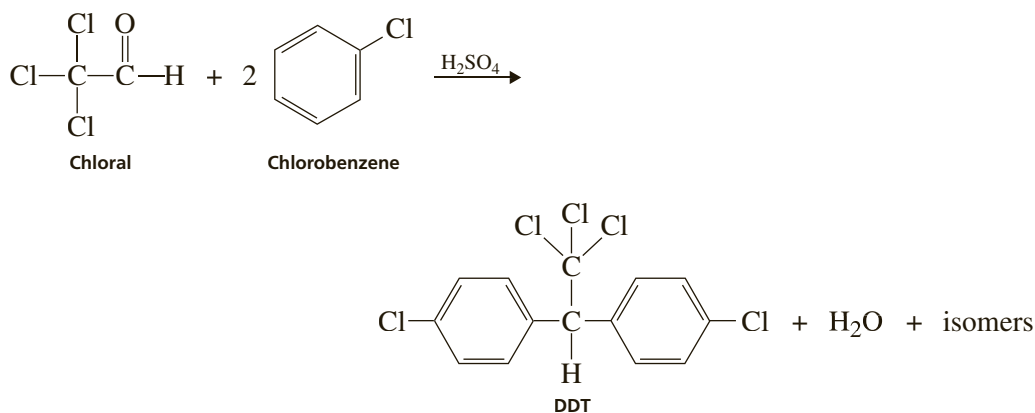


At the time DDT was introduced, it was an important boon to humanity. It was effective in controlling lice, fleas, and malaria-carrying mosquitoes and thus helped control human and animal disease. The use of DDT rapidly spread to the control of hundreds of insects that damage fruit, vegetable, and grain crops.

Pesticides that persist in the environment for a long time after application are called **hard pesticides**. Beginning in the 1960s, some of the harmful effects of such hard pesticides as DDT and the other chlorocarbon materials became known. DDT is a fat-soluble material and is therefore likely to collect in the fat, nerve, and brain tissues of animals. The concentration of DDT in tissues increases in animals high in the food chain. Thus, birds that eat poisoned insects accumulate large quantities of

DDT. Animals that feed on the birds accumulate even more DDT. In birds, at least two undesirable effects of DDT have been recognized. First, birds whose tissues contain large amounts of DDT have been observed to lay eggs having shells too thin to survive until young birds are ready to hatch. Second, large quantities of DDT in the tissues seem to interfere with normal reproductive cycles. The massive destruction of bird populations that sometimes occurred after heavy spraying with DDT became an issue of great concern. The brown pelican and the bald eagle were placed in danger of extinction. The use of chlorocarbon insecticides was identified as the principal reason for the decline in the numbers of these birds.

Because DDT is chemically inert, it persists in the environment without decomposing to harmless materials. It can decompose very slowly, but the decomposition products are every bit as harmful as DDT itself. Consequently, each application of DDT means that still more DDT will pass from species to species—from food source to predator—until it concentrates in the higher animals, possibly endangering their existence. Even humans may be threatened. As a result of evidence of the harmful effects of DDT, the Environmental Protection Agency (EPA) banned general use of DDT in the early 1970s; it may still be used for certain purposes, although permission of the EPA is required. In 1974, the EPA granted permission to use DDT against the tussock moth in the forests of Washington and Oregon.


Because the life cycles of insects are short, they can evolve an immunity to insecticides within a short period. As early as 1948, several strains of DDT-resistant insects were identified. Today, the malaria-bearing mosquitoes are almost completely resistant to DDT, an ironic development. Other chlorocarbon insecticides were developed to use as alternatives to DDT against resistant insects. Examples of these chlorocarbon materials include Dieldrin, Aldrin, Chlordane, and the substances whose structures are shown here. Heptachlor and Mirex are prepared using Diels–Alder reactions.

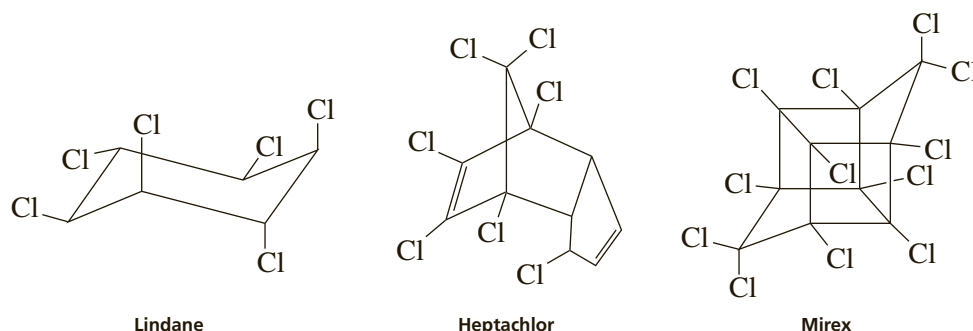

In spite of structural similarity, Chlordane and Heptachlor behave differently than DDT, Dieldrin, and Aldrin. Chlordane, for instance, is short-lived and less toxic to mammals. Nevertheless, all the chlorocarbon insecticides have been the objects of much suspicion. A ban on the use of Dieldrin and Aldrin has also been ordered by the EPA. In addition, strains of insects resistant to Dieldrin, Aldrin, and other materials have been observed. Some insects become addicted to a chlorocarbon insecticide and thrive on it!


The problems associated with chlorocarbon materials have led to the development of “soft” insecticides. These usually are organophosphorus or carbamate derivatives, and they are characterized by a short lifetime before they are decomposed to harmless materials in the environment.

The organic structures of some organophosphorus insecticides are shown here.

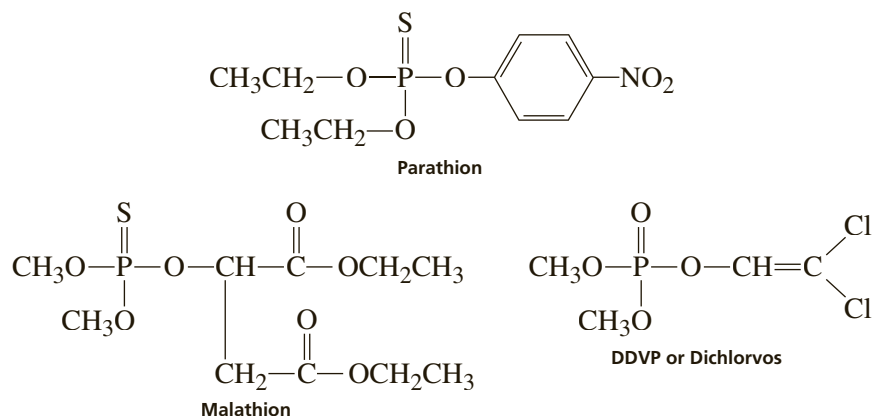

Parathion and Malathion are used widely for agriculture. DDVP is contained in “pest strips,” which are used to combat household insect pests. The organophosphorus materials do not persist in the environment, so they are not passed between species up the food chain, as the chlorocarbon compounds are. However, the organophosphorus compounds are highly toxic to humans. Some migrant and other agricultural workers have lost their lives because of accidents involving these materials. Stringent safety precautions must be applied when organophosphorus insecticides are being used.

The carbamate derivatives, including Carbaryl, tend to be less toxic than the organophosphorus compounds. They are also readily degraded to harmless materials. Nevertheless, insects resistant to soft insecticides have also been observed. Furthermore, the organophosphorus and carbamate derivatives destroy many more nontarget pests than the chlorocarbon compounds do. The danger to earthworms, mammals, and birds is very high.

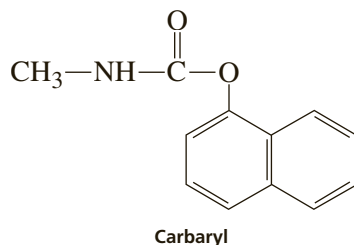


ALTERNATIVES TO INSECTICIDES


Several alternatives to the massive application of insecticides have recently been explored. Insect attractants, including the pheromones (see the essay preceding Experiment 45), have been used in localized traps. Such methods have been effective against the gypsy moth. A “confusion technique,” whereby a pheromone is sprayed into the air in such high concentrations that male insects are no longer able to locate females, has been studied. These methods are specific to the target pest and do not cause repercussions in the general environment.

Recent research has focused on using an insect’s own biochemical processes to control pests. Experiments with **juvenile hormone** have shown promise. Juvenile hormone is one of three internal secretions used by insects to regulate growth and metamorphosis from larva to pupa and thence to the adult. At certain stages in the metamorphosis from larva to pupa, juvenile hormone must be secreted; at other stages it must be absent, or the insect will either develop abnormally or fail to mature. Juvenile hormone is important in maintaining the juvenile, or larval, stage of the growing insect. The male cecropia moth, which is the mature form of the

silkworm, has been used as a source of juvenile hormone. The structure of the cecropia juvenile hormone is shown below. This material has been found to prevent the maturation of yellow-fever mosquitoes and human body lice. Because insects are not expected to develop a resistance to their own hormones, it is hoped that insects will be unlikely to develop a resistance to juvenile hormone.

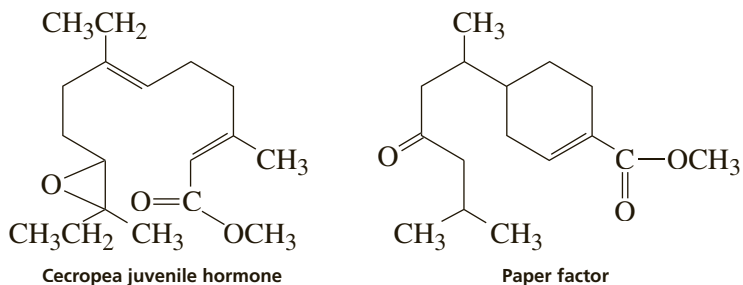

Although it is very difficult to get enough of the natural substance for use in agriculture, synthetic analogues have been prepared, that have been shown to be similar in properties and effectiveness to the natural substance. A substance has been found in the American balsam fir (*Abies balsamea*) known as **paper factor**. Paper factor is active against the linden bug, *Pyrrhocoris apterus*, a European cotton pest. This substance is merely one of thousands of terpenoid materials synthesized by the fir tree. Other terpenoid substances are being investigated as potential juvenile hormone analogues.

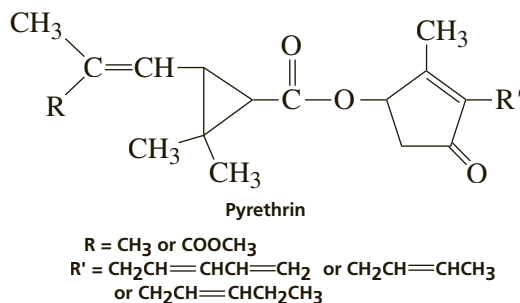

Certain plants are capable of synthesizing substances that protect them against insects. Included among these natural insecticides are the **pyrethrins** and derivatives of **nicotine**.

The search for environmentally suitable means of controlling agricultural pests continues with a great sense of urgency. Insects cause billions of dollars of damage to food crops each year. With food becoming increasingly scarce and with the world's population growing at an exponential rate, preventing such losses to food crops is absolutely essential.

REFERENCES

- Berkoff, C. E. Insect Hormones and Insect Control. *J. Chem. Educ.* **1971**, *48*, 577.
- Bowers, W. S; Nishida, R. Juvocimenes: Potent Juvenile Hormone Mimics from Sweet Basil. *Science* **1980**, *209*, 1030.
- Carson, R. *Silent Spring*; Houghton Mifflin: Boston. 1962.


- Keller, E. The DDT Story. *Chemistry* **1970**, 43 (Feb), 8.
- O'Brien, R. D. *Insecticides: Action and Metabolism*; Academic Press: New York, 1967.
- Peakall, D. B. Pesticides and the Reproduction of Birds. *Sci Am.* **1970**, 222 (Apr), 72.
- Saunders, H. J. New Weapons against Insects. *Chem. Eng. News* **1975**, 53 (Jul 28), 18.
- Williams, C. M. Third-Generation Pesticides. *Sci. Am.* **1967**, 217 (Jul), 13.
- Williams, W. G.; Kennedy, G. G.; Yamamoto, R. T.; Thacker, J. D.; Bordner, J. 2-Tridecanone: A Naturally Occurring Insecticide from the Wild Tomato. *Science* **1980**, 207, 888.

The Diels–Alder Reaction of Cyclopentadiene with Maleic Anhydride


Diels–Alder reaction

Fractional distillation

Cyclopentadiene and maleic anhydride react readily in a Diels–Alder reaction to form the adduct, *cis*-norbornene-5,6-*endo*-dicarboxylic anhydride:

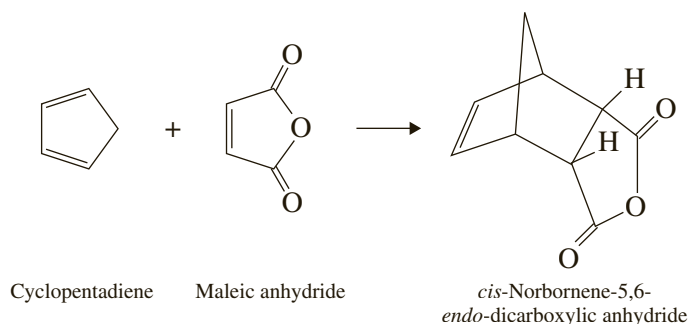

Because two molecules of cyclopentadiene can also undergo a Diels–Alder reaction to form dicyclopentadiene, it is not possible to store cyclopentadiene in the monomeric form. Therefore, it is necessary to first “crack” dicyclopentadiene to produce cyclopentadiene for use in this experiment. This will be accomplished by heating the dicyclopentadiene to a boil and collecting the cyclopentadiene as it is formed by fractional distillation. To keep it from dimerizing, the cyclopentadiene must be kept cold and used fairly soon.

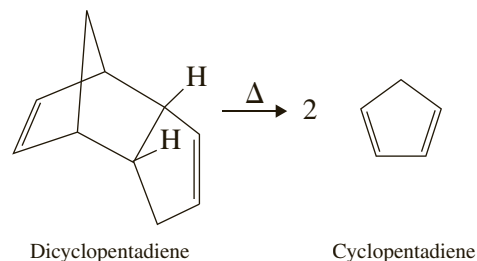

REQUIRED READING

Review: Technique 11 Crystallization, Section 11.4

New: Essay Diels–Alder Reaction and Insecticides

SPECIAL INSTRUCTIONS

The cracking of dicyclopentadiene should be performed by the instructor or laboratory assistant. If a flame is used for this, be sure that there are no leaks in the system, because both cyclopentadiene and the dimer are highly flammable. The procedure provides enough cyclopentadiene for about 50 students.

SUGGESTED WASTE DISPOSAL

Dispose of the mother liquor from the crystallization in the container designated for nonhalogenated organic solvents.

NOTES TO THE INSTRUCTOR

Working in a hood, assemble a fractional distillation apparatus, as shown in the figure. Glassware with a joint size of $\text{F } 19/22$ or larger should be used. If smaller glassware is used, the fractionating column may not be long enough to achieve the necessary separation. Although the required temperature control can best be obtained with a microburner, using a heating mantle, aluminum block, or sand bath lessens the possibility of a fire occurring. Place several boiling stones and 15 mL of dicyclopentadiene in the 50-mL distilling flask. Control the heat source so that the cyclopentadiene distills at 40–43°C. (If a sand bath is used, the temperature should be 190–200°C, and it may be necessary to cover the sand bath and distilling flask with aluminum foil.) After 30–45 minutes, 6–7 mL of cyclopentadiene should be collected, and the distillation can be stopped. If the cyclopentadiene is cloudy, dry the liquid over granular anhydrous sodium sulfate. Store the product in a sealed container and keep it cooled in an ice-water bath until all students have taken their portions.


PROCEDURE

Preparation of the Adduct

To a Craig tube add 0.100 g of maleic anhydride and 0.40 mL of ethyl acetate. Without inserting the plug, shake the tube gently to dissolve the solid (slight heating in a warm water bath may be necessary). Add 0.40 mL of ligroin (bp 60–90°C) and shake the tube gently to mix the solvents and reactant thoroughly. Add 0.10 mL of cyclopentadiene and mix thoroughly by shaking until no visible layers of liquid are present. Because this reaction is exothermic, the temperature of the mixture will likely become high enough to keep the product in solution. However, if a solid does form at this point, it will be necessary to heat the mixture gently in a warm water bath to dissolve any solids present. If necessary, add a drop of ethyl acetate to help dissolve the solid, and, again, heat the mixture gently.

Crystallization of Product

Allow the mixture to cool slowly to room temperature by placing the Craig tube in a 10-mL Erlenmeyer flask that has been filled with about 8 mL of water at 50–60°C. The inner plug of the Craig tube should be inserted to prevent

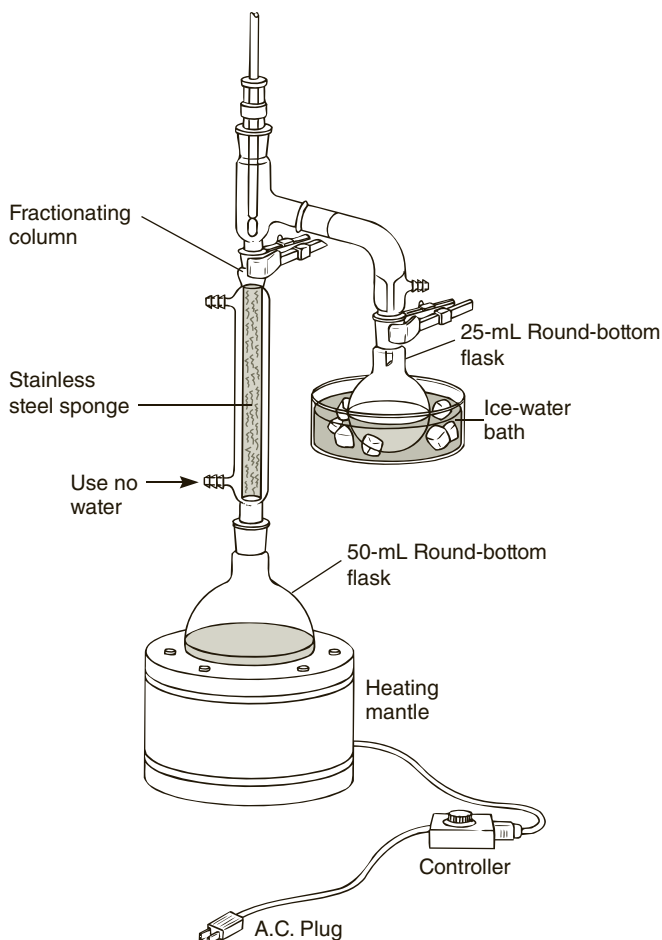

Fractional distillation apparatus for cracking dicyclopentadiene.


evaporation of the solvent. Better crystal formation can be achieved by seeding the solution before it cools to room temperature. To seed the solution, dip a spatula or glass stirring rod into the solution after it has cooled for about 5 minutes. Allow the solvent to evaporate so that a small amount of solid forms on the surface of the spatula or glass rod. Place the spatula or stirring rod back into the solution for a few seconds to induce crystallization. When crystallization is complete at room temperature, cool the mixture in an ice bath for several minutes.

Isolate the crystals from the Craig tube by centrifugation (see Technique 8, Section 8.7, and Figure 8.11) and allow the crystals to air-dry. Determine the weight and the melting point (164°C).

At the option of the instructor, obtain the infrared spectrum using the dry film method (Technique 25, Section 25.4) or by an alternative method (Technique 25, Section 25.5). Compare your infrared spectrum with the one reproduced here. Calculate the percentage yield and submit the product to the instructor in a labeled vial.


MOLECULAR MODELING (OPTIONAL)

In the reaction of cyclopentadiene with maleic anhydride, two products are possible: the *endo* product and the *exo* product.

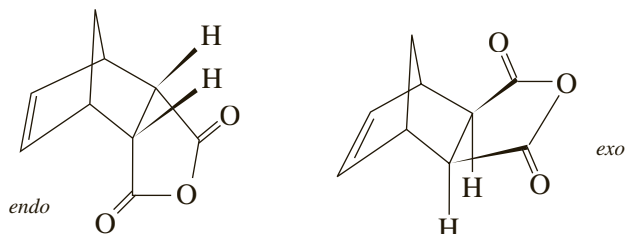


Calculate the heats of formation for both of these products to determine which is the expected **thermodynamic product** (product of lowest energy). Perform the calculations at the AM1 level with a geometry optimization. The actual product of the Diels–Alder reaction is the *endo* product; is this the thermodynamic product? Display a space-filling model for each structure. Which one appears most crowded?

Woodward and Hoffmann have pointed out that the diene is the electron donor and the dienophile the electron acceptor in this reaction. In accordance with this idea, dienes that have electron-donating groups are more reactive than those

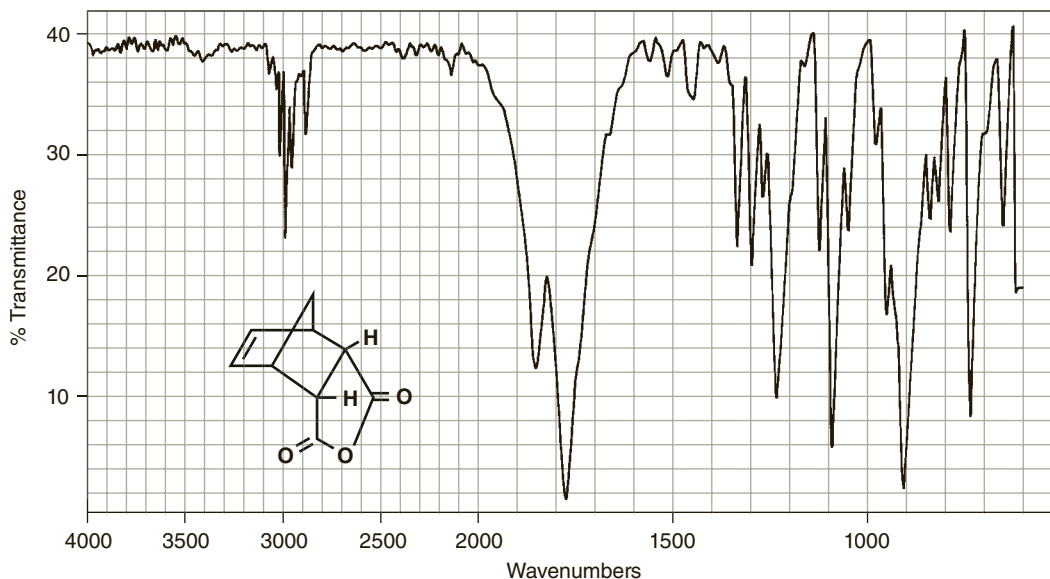

Infrared spectrum of *cis*-norbornene-5,6-*endo*-dicarboxylic anhydride, KBr.

without, and dienophiles with electron-withdrawing groups are most reactive. Using the reasoning of frontier molecular orbital theory (see the essay “Computational Chemistry”), the electrons in the HOMO of the diene will be placed into the LUMO of the dienophile when reaction occurs. Using the AM1 level, calculate the HOMO surface for the diene (cyclopentadiene) and the LUMO surface for the dienophile (maleic anhydride). Display the two simultaneously on the screen in the orientations that will lead to the *endo* and *exo* products.

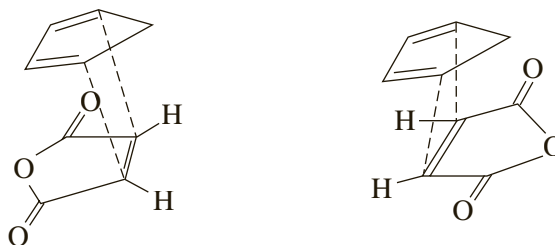

Woodward and Hoffmann suggested that the orientation that leads to the largest degree of constructive overlap between the two orbitals (HOMO and LUMO) is the orientation that would lead to the product. Do you agree?


Depending on the capability of your software, it may be possible to determine the geometrics (and energies) of the transition states that lead to each product. Your instructor will have to show you how to do this.

QUESTIONS

1. Draw a structure for the *exo* product formed by cyclopentadiene and maleic anhydride.
2. Because the *exo* form is more stable than the *endo* form, why is the *endo* product formed almost exclusively in this reaction?
3. In addition to the main product, what are two side reactions that could occur in this experiment?
4. The infrared spectrum of the adduct is given in this experiment. Interpret the principal peaks.

The Diels–Alder Reaction with Anthracene-9-methanol

Green Chemistry
Diels–Alder reaction
Hydrophobic effect
Spectroscopy

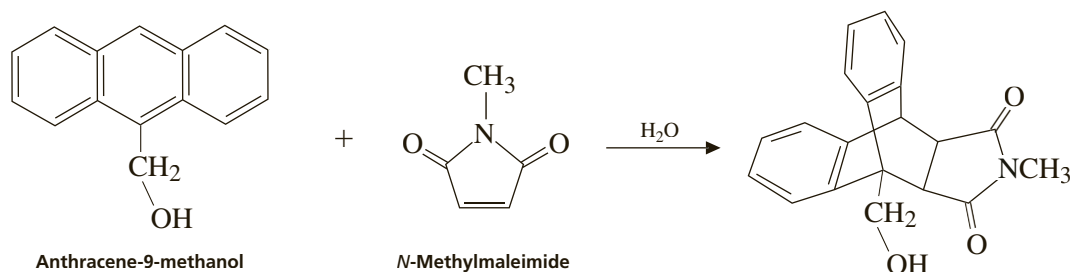

This experiment demonstrates Green Chemistry through the Diels–Alder reaction, which is an important reaction in organic chemistry because it is an important method of ring formation. The “green” components of this experiment include attention to atom economy and waste reduction, but the most important “green” aspect is the use of water as the solvent. Not only is water an environmentally benign solvent, but it also actually improves other aspects of this reaction due to hydrophobic solvent effects.

The *hydrophobic effect* is the property that nonpolar molecules tend to self-associate in the presence of aqueous solution. Two explanations have been advanced to explain why the hydrophobic effect increases the rate of reaction for selected Diels–Alder reactions. The first is that the activated complex is somewhat polar; it is stabilized by hydrogen-bonding, which makes the reaction go faster. The second is that the hydrophobic effect acts to force the two reagents together with a solvation shell and to increase the interaction between them.

SUGGESTED WASTE DISPOSAL

All aqueous waste can be disposed of in a waste container designated for nonhalogenated aqueous waste.

SAFETY PRECAUTIONS

N-Methylmaleimide is corrosive and should be handled with care. Gloves should be worn.

PROCEDURE

Reaction

In a 50-mL round-bottom flask equipped with a stir bar, add 0.066 g of anthracene-9-methanol. Using a 25-mL graduated cylinder, add 25 mL of de-ionized water. Note that anthracene-9-methanol is insoluble in water. Add 0.070 g of *N*-methylmaleimide to the mixture, and fit the flask to a water-cooled condenser. Heat the mixture until it is boiling under reflux, and allow the reaction to continue boiling for 90 minutes while stirring.

Isolation


Remove the heat, and allow the reaction to cool to room temperature (without removing the condenser). Chill the flask in an ice bath for 5 minutes, and collect the precipitate by vacuum filtration using a Hirsch funnel. Allow the solid to dry in the Hirsch funnel, under vacuum, for 15 minutes. Collect crystals on a watch glass, and allow them to dry overnight.

Analysis and Report

Determine the weight of your product, and obtain the melting point range (literature value = 232–235°C). Determine the proton and carbon nuclear magnetic resonance spectra of the product. Include the NMR spectra with your report, along with an interpretation of the peaks and splitting patterns. Be sure to also include your weight percentage recovery calculation. Submit your sample in a properly labeled vial.


NMR SPECTRAL ANALYSIS

If you examine the complete structural formula of the Diels–Alder product in this experiment, you will see that the hydrogen atoms labeled **f** and **g** are not equivalent! Owing to the position of the ring labeled **C** in the formula (it is shown tilted to the left), the overall ring system is not symmetrical. The ring labeled **A** is not equivalent to ring **B**.

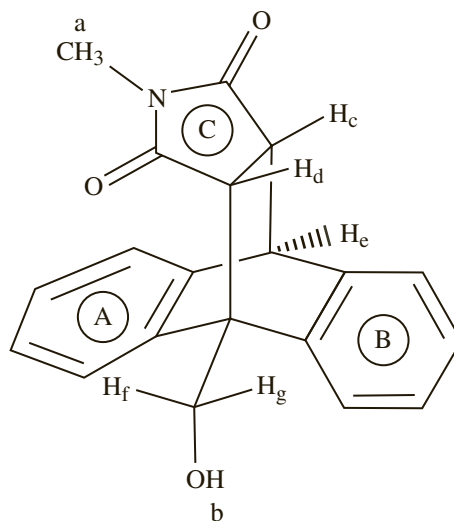

An examination of the Newman projection of the $-\text{CH}_2\text{OH}$ part of the molecule shows clearly that there is no conformation in which proton **f** becomes equivalent to

proton **g**. The two protons are thus **diastereotopic** (see Technique 26, Section 26.16 for a discussion of diastereotopic protons). The two protons appear at different values of chemical shift.

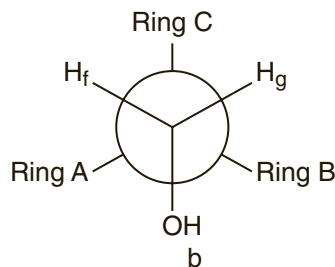

Since proton **f** and proton **g** are not equivalent, they will couple to each other in the NMR spectrum. This coupling can be seen in the expanded NMR spectrum (Section 26.16). You will also be able to see coupling between the –OH proton (**b**) and protons **f** and **g**. Because proton **f** is not equivalent to proton **g**, J_{bf} may not be equal to J_{bg} .

REFERENCES

- Engberts, J. B. F. N. Diels–Alder Reactions in Water: Enforced Hydrophobic Interaction and Hydrogen Bonding. *Pure Appl. Chem.* **1995**, *67*, 823–828.
- Rideout, D. C.; Breslow, R. Hydrophobic Acceleration of Diels–Alder Reactions. *J. Am. Chem. Soc.* **1980**, *102*, 7817–7818.

Photoreduction of Benzophenone and Rearrangement of Benzpinacol to Benzopinacolone

Photochemistry


Photoreduction

Energy transfer


Pinacol rearrangement

This experiment consists of two parts. In the first part (Experiment 55A), benzophenone will be subjected to **photoreduction**, a dimerization brought about by exposing a solution of benzophenone in isopropyl alcohol to natural sunlight. The product of this photoreaction is benzpinacol. In the second part (Experiment 55B), benzpinacol will be induced to undergo an acid-catalyzed rearrangement called the **pinacol rearrangement**. The product of the rearrangement is benzopinacolone.

Experiment 55A

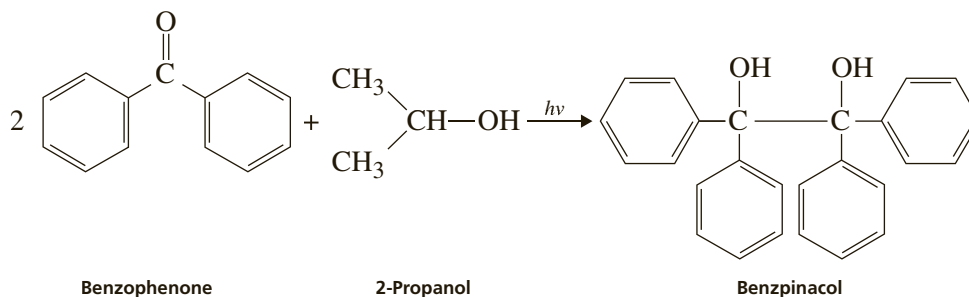


Experiment 55B

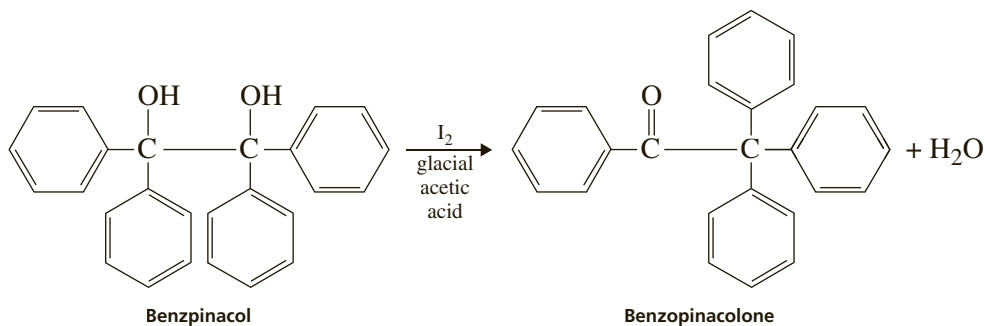


55A EXPERIMENT 55A

Photoreduction of Benzophenone

The photoreduction of benzophenone is one of the oldest and most thoroughly studied photochemical reactions. Early in the history of photochemistry, it was discovered that solutions of benzophenone are unstable in light when certain solvents are used. If benzophenone is dissolved in a “hydrogen-donor” solvent, such as 2-propanol, and exposed to ultraviolet light, $h\nu$, an insoluble dimeric product, benzpinacol, will form.

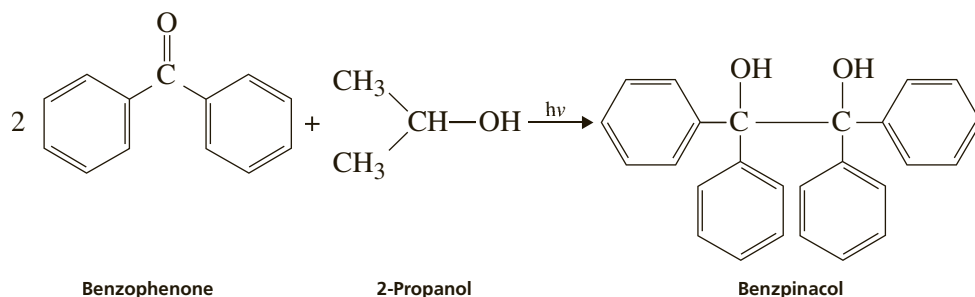

To understand this reaction, let's review some simple photochemistry as it relates to aromatic ketones. In the typical organic molecule, all the electrons are paired in the occupied orbitals. When such a molecule absorbs ultraviolet light of the appropriate wavelength, an electron from one of the occupied orbitals, usually the one of highest energy, is excited to an unoccupied molecular orbital, usually to the one of lowest energy. During this transition, the electron must retain its spin value, because during an electronic transition a change of spin is forbidden by the laws of quantum mechanics. Therefore, just as the two electrons in the highest occupied orbital of the molecule originally had their spins paired (opposite), so they will retain paired spins in the first electronically excited state of the molecule. This is true even though the two electrons will be in *different* orbitals after the transition.

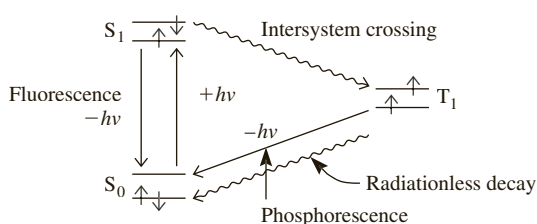

Electronic states of a typical molecule and the possible interconversions. In each state (S_0 , S_1 , T_1), the lower line represents the highest occupied orbital, and the upper line represents the lowest unoccupied orbital of the unexcited molecule. Straight lines represent processes in which a photon is absorbed or emitted. Wavy lines represent radiationless processes—those that occur without emission or absorption of a photon.


This first excited state of a molecule is called a **singlet state** (S_1) because its spin multiplicity ($2S + 1$) is 1. The original unexcited state of the molecule is also a singlet state because its electrons are paired, and it is called the **ground-state singlet state** (S_0) of the molecule.

The excited state singlet S_1 may return to the ground state S_0 by reemission of the absorbed photon of energy. This process is called **fluorescence**. Alternatively, the excited electron may undergo a change of spin to give a state of higher multiplicity, the excited **triplet state**, so called because its spin multiplicity ($2S + 1$) is 3. The conversion from the first excited singlet state to the triplet state is called **intersystem crossing**. Because the triplet state has a higher multiplicity, it inevitably has a lower energy state than the excited singlet state (Hund's Rule). Normally, this change of spin (intersystem crossing) is a process forbidden by quantum mechanics, just as a direct excitation of the ground state (S_0) to the triplet state (T_1) is forbidden. However, in


those molecules in which the singlet and triplet states lie close to one another in energy, the two states inevitably have several overlapping vibrational states—that is, states in common—a situation that allows the “forbidden” transition to occur. In many molecules in which S_1 and T_1 have similar energy ($\Delta E < 10$ Kcal/mole), intersystem crossing occurs faster than fluorescence, and the molecule is rapidly converted from its excited singlet state to its triplet state. In benzophenone, S_1 undergoes intersystem crossing to T_1 with a rate of $k_{isc} = 10^{10} \text{ sec}^{-1}$, meaning that the lifetime of S_1 is only 10^{-10} second. The rate of fluorescence for benzophenone is $k_f = 10^6 \text{ sec}^{-1}$, meaning that intersystem crossing occurs at a rate that is 10^4 times faster than fluorescence. Thus, the conversion of S_1 to T_1 in benzophenone is essentially a quantitative process. In molecules that have a wide energy gap between S_1 and T_1 , this situation would be reversed. As you will see shortly, the naphthalene molecule presents a reversed situation.

Because the excited triplet state is lower in energy than the excited singlet state, the molecule cannot easily return to the excited singlet state. Nor can it easily return to the ground state by returning the excited electron to its original orbital. Once again, the transition $T_1 \rightarrow S_0$ would require a change of spin for the electron, and this is a forbidden process. Hence, the triplet excited state usually has a long lifetime (relative to other excited states) because it generally has nowhere to which it can easily go. Even though the process is forbidden, the triplet T_1 may eventually return to the ground state (S_0) by a process called a **radiationless transition**. In this process, the excess energy of the triplet is lost to the surrounding solution as heat, thereby “relaxing” the triplet back to the ground state (S_0). This process is the study of much current research and is not well understood. In the second process, in which a triplet state may revert to the ground state, **phosphorescence**, the excited triplet emits a photon to dissipate the excess energy and returns directly to the ground state. Although this process is “forbidden,” it nevertheless occurs when there is no other open pathway by which the molecule can dissipate its excess energy. In benzophenone, radiationless decay is the faster process, with rate $k_d = 10^5 \text{ sec}^{-1}$, and phosphorescence, which is not observed, has a lower rate of $k_p = 10^2 \text{ sec}^{-1}$.

Benzophenone is a ketone. Ketones have *two* possible excited singlet states and, consequently, two excited triplet states as well. This occurs because two relatively low-energy transitions are possible in benzophenone. It is possible to excite one of the π electrons in the carbonyl π bond to the lowest-energy unoccupied orbital, a π^* orbital. It is also possible to excite one of the unbonded or n electrons on oxygen to the same orbital. The first type of transition is called a $\pi-\pi^*$ transition, whereas the second is called an $n-\pi^*$ transition. In the figure, these transitions and the states that result are illustrated pictorially.

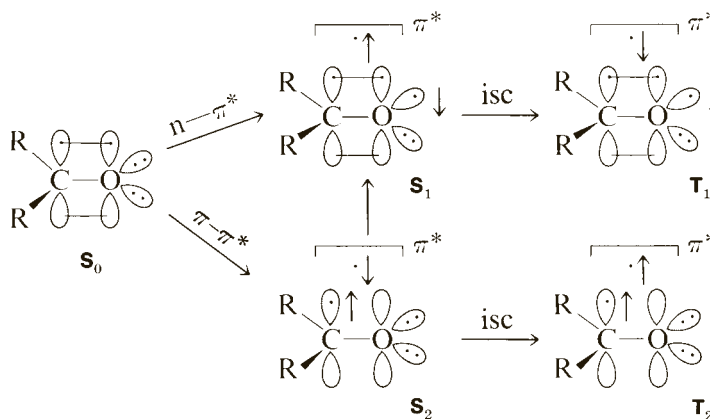




$n-\pi^*$ and $\pi-\pi^*$ transitions for ketones.

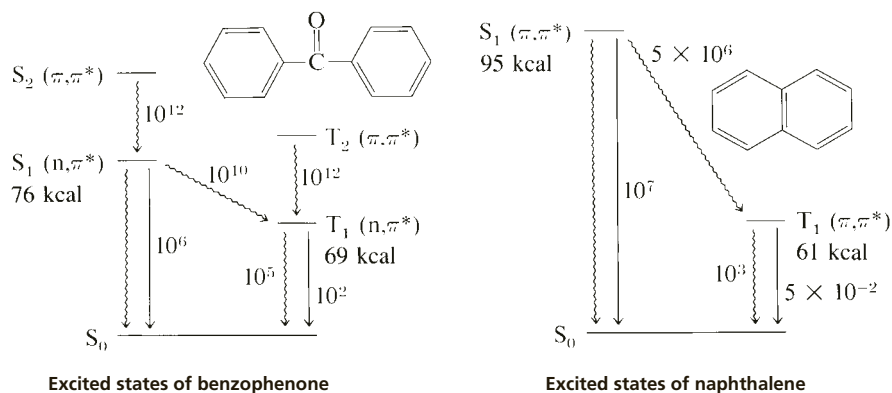

Excited energy states of benzophenone and naphthalene.

Spectroscopic studies show that for benzophenone and most other ketones, the $n-\pi^*$ excited states S_1 and T_1 are of lower energy than the $\pi-\pi^*$ excited states. An energy diagram depicting the excited states of benzophenone (along with one that depicts those of naphthalene) is shown.

It is now known that the photoreduction of benzophenone is a reaction of the $n-\pi^*$ triplet state (T_1) of benzophenone. The $n-\pi^*$ excited states have radical character at the carbonyl oxygen atom because of the unpaired electron in the nonbonding orbital. Thus, the radical-like and energetic T_1 excited-state species can abstract a hydrogen atom from a suitable donor molecule to form the diphenylhydroxymethyl radical. Two of these radicals, once formed, may couple to form benzpinacol. The complete mechanism for photoreduction is outlined in the steps that follow.

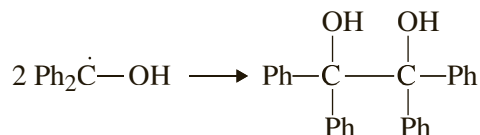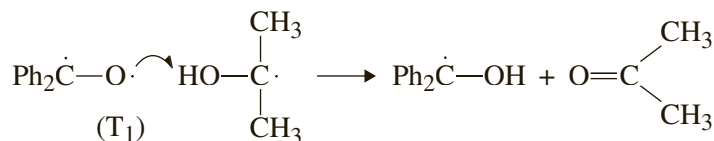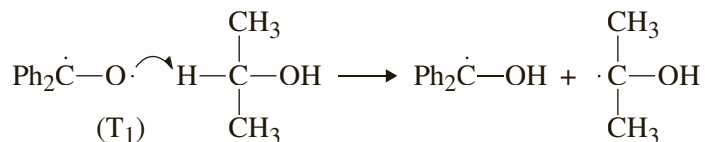

Many photochemical reactions must be carried out in a quartz apparatus because they require ultraviolet radiation of shorter wavelengths (higher energy) than the wavelengths that can pass through Pyrex. Benzophenone, however, requires radiation of approximately 350 nm to become excited to its $n-\pi^*$ singlet state S_1 , a wavelength that readily passes through Pyrex. In the figure shown below, the ultraviolet absorption spectra of benzophenone and naphthalene are given. Superimposed on their spectra are two curves, which show the wavelengths that can be transmitted by Pyrex and quartz, respectively. Pyrex will not allow any radiation of wavelength shorter than approximately 300 nm to pass, whereas quartz allows wavelengths as



short as 200 nm to pass. Thus, when benzophenone is placed in a Pyrex flask, the only electronic transition possible is the $n-\pi^*$ transition, which occurs at 350 nm.

However, even if it were possible to supply benzophenone with radiation of the appropriate wavelength to produce the second excited singlet state of the molecule, this singlet would rapidly convert to the lowest singlet state (S_1). The state S_2 has a lifetime of less than 10^{-12} second. The conversion process $S_2 \rightarrow S_1$ is called an **internal conversion**. Internal conversions are processes of conversion between excited states of the same multiplicity (singlet–singlet or triplet–triplet), and they usually are rapid. Thus, when an S_2 or T_2 is formed, it readily converts to S_1 or T_1 , respectively. As a consequence of their short lifetimes, little is known about the properties or the exact energies of S_2 and T_2 of benzophenone.


Energy Transfer

Using a simple **energy-transfer** experiment, one can show that the photoreduction of benzophenone proceeds via the T_1 excited state of benzophenone rather than the S_1 excited state. If naphthalene is added to the reaction, the photoreduction is stopped because the excitation energy of the benzophenone triplet is transferred to naphthalene. The naphthalene is said to have **quenched** the reaction. This occurs in the following way.

When the excited states of molecules have long enough lifetimes, they often can transfer their excitation energy to another molecule. The mechanisms of these transfers are complex and cannot be explained here; however, the essential requirements can be outlined. First, for two molecules to exchange their respective states of excitation, the process must occur with an overall decrease in energy. Second, the spin multiplicity of the total system must not change. These two features can be illustrated by the two most common examples of energy transfer: singlet transfer and triplet transfer. In these two examples, the superscript 1 denotes an excited singlet state, the superscript 3 denotes a triplet state, and the subscript 0 denotes a ground-state molecule. The designations A and B represent different molecules.

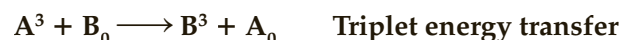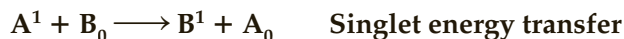

In singlet energy transfer, excitation energy is transferred from the excited singlet state of A to a ground-state molecule of B, converting B to its excited singlet state and returning A to its ground state. In triplet energy transfer, there is a similar interconversion of excited state and ground state. Singlet energy is transferred

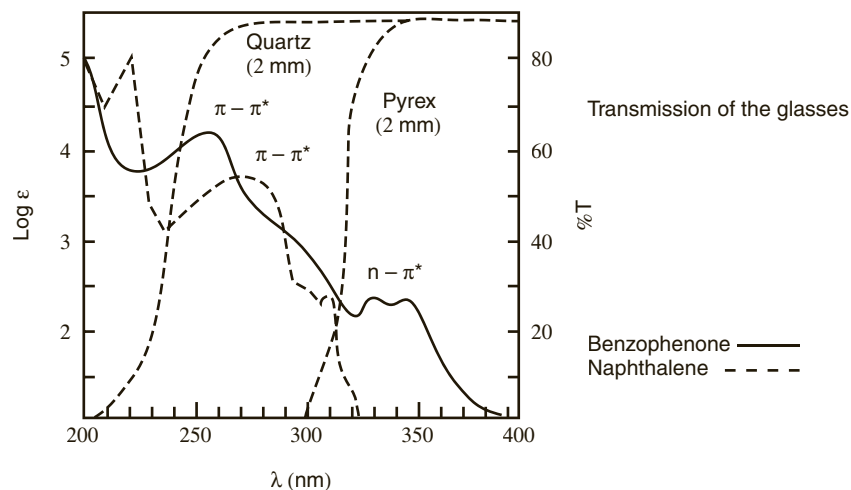

Ultraviolet absorption spectra for benzophenone and naphthalene.

through space by a dipole–dipole coupling mechanism, but triplet energy transfer requires the two molecules involved in the transfer to collide. In the usual organic medium, about 10^9 collisions occur per second. Thus, if a triplet state A^3 has a lifetime longer than 10^{-9} second, and if an acceptor molecule B_0 , which has a lower triplet energy than that of A^3 , is available, energy transfer can be expected. If the triplet A^3 undergoes a reaction (such as photoreduction) at a rate lower than the rate of collisions in the solution, and if an acceptor molecule is added to the solution, the reaction can be *quenched*. The acceptor molecule, which is called a **quencher**, deactivates, or “quenches,” the triplet before it has a chance to react. Naphthalene has the ability to quench benzophenone triplets in this way and to stop the photoreduction.

Naphthalene cannot quench the excited-state singlet S_1 of benzophenone because its own singlet has an energy (95 kcal/mol) that is higher than the energy of benzophenone (76 kcal/mol). In addition, the conversion $S_1 \rightarrow T_1$ is rapid (10^{-10} second) in benzophenone. Thus, naphthalene can intercept only the triplet state of benzophenone. The triplet excitation energy of benzophenone (69 kcal/mol) is transferred to naphthalene ($T_1 = 61$ kcal/mol) in an exothermic collision. Finally, the naphthalene molecule does not absorb light of the wavelengths transmitted by Pyrex (see the spectra shown here); therefore, benzophenone is not inhibited from absorbing energy when naphthalene is present in solution. Thus, because naphthalene quenches the photoreduction reaction of benzophenone, we can infer that this reaction proceeds via the triplet state T_1 of benzophenone. If naphthalene did not quench the reaction, the singlet state of benzophenone would be indicated as the reactive intermediate. In the following experiment, the photoreduction of benzophenone is attempted both in the presence and in the absence of added naphthalene.

REQUIRED READING

Review: Technique 8 Filtration, Section 8.3

SPECIAL INSTRUCTIONS

This experiment may be performed concurrently with some other experiment. It requires only 15 minutes during the first laboratory period and only about 15 minutes in a subsequent laboratory period about 1 week later (or at the end of the laboratory period if you use a sunlamp).

Using Direct Sunlight

It is important that the reaction mixture be left where it will receive direct sunlight. If it does not, the reaction will be slow and may need more than 1 week for completion. It is also important that the room temperature not be too low or the benzophenone will precipitate. If you perform this experiment in the winter and the laboratory is not heated at night, you must shake the solutions every morning to redissolve the benzophenone. Benzopinacol should not redissolve easily.

Using a Sunlamp

If you wish, you may use a 275-W sunlamp instead of direct sunlight. Place the lamp in a hood that has had its window covered with aluminum foil (shiny side in). The lamp (or lamps) should be mounted in a ceramic socket attached to a ring stand with a three-pronged clamp.

CAUTION

The purpose of the aluminum foil is to protect the eyes of people in the laboratory. You should not view a sunlamp directly or damage to the eyes may result. Take all possible viewing precautions.

Attach samples to a ring stand placed at least 18 inches from the sunlamp. Placing them at this distance will avoid their being heated by the lamp. Heating may cause loss of the solvent. It is a good idea to agitate the samples every 30 minutes. With a sunlamp, the reaction will be complete in 3–4 hours.

SUGGESTED WASTE DISPOSAL

Dispose of the filtrate from the vacuum filtration procedure in the container designated for nonhalogenated organic wastes.

PROCEDURE

Label two 13 × 100-mm test tubes near the top of the tubes. The labels should have your name and “No. 1” and “No. 2” written on them. Place 0.50 g of benzophenone in the first tube. Place 0.50 g of benzophenone and 0.05 g of naphthalene in the second tube. Add about 2 mL of 2-propanol (isopropyl alcohol) to each tube and warm them in a beaker of warm water to dissolve the solids. When the solids have dissolved, add one small drop (Pasteur pipette) of glacial acetic acid to each tube and then fill each tube nearly to the top with more 2-propanol. Stopper the tubes tightly with rubber stoppers, shake them well, and place them in a beaker on a windowsill where they will receive direct sunlight.

NOTE: You may be directed by your instructor to use a sunlamp instead of direct sunlight (see “Special Instructions”).

The reaction requires about 1 week for completion (3 hours with a sunlamp). If the reaction has occurred during this period, the product will have crystallized from the solution. Observe the result in each test tube. Collect the product by vacuum filtration using a small Büchner or Hirsch funnel (Technique 8, Section 8.3) and allow it to dry. Weigh the product and determine its melting point and percentage yield. At the option of the instructor, obtain the infrared spectrum using the dry film method (Technique 25, Section 25.4) or by an alternative method (Technique 25, Section 25.5). Submit the product to the instructor in a labeled vial along with the report.


REFERENCE

Vogler, A.; Kunkely, H. Photochemistry and Beer. *J. Chem. Educ.*, 1982, 59 (January), 25.

55B EXPERIMENT 55B

Synthesis of β -Benzopinacolone: The Acid-Catalyzed Rearrangement of Benzpinacol

The ability of carbocations to rearrange represents an important concept in organic chemistry. In this experiment, the benzpinacol, prepared in Experiment 55A, will rearrange to **benzopinacolone (2,2,2-triphenylacetophenone)** under the influence of iodine in glacial acetic acid.

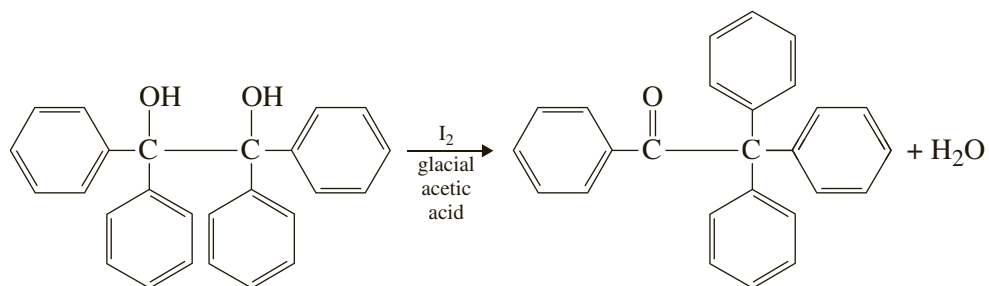

The product is isolated as a crystalline white solid. Benzopinacolone is known to crystallize in two crystalline forms, each with a different melting point. The **alpha** form has a melting point of 206–207°C, whereas the **beta** form melts at 182°C. The product formed in this experiment is the β -benzopinacolone.

REQUIRED READING

- Review:* Technique 7 Reaction Methods, Section 7.2
Technique 11 Crystallization: Purification of Solids, Section 11.3
Technique 25 Infrared Spectroscopy, Part B
Technique 26 Nuclear Magnetic Resonance Spectroscopy, Part B

Before beginning this experiment, you should read the material dealing with carbocation rearrangements in your lecture textbook.

SPECIAL INSTRUCTIONS

This experiment requires little time and can be coscheduled with another short experiment.

SUGGESTED WASTE DISPOSAL

All organic residues must be placed in the appropriate container for nonhalogenated organic waste.

PROCEDURE


In a 25-mL round-bottom flask, add 5 mL of a 0.015 M solution of iodine dissolved in glacial acetic acid. Add 1 g of benzpinacol and attach a water-cooled condenser. Using a small heating mantle, allow the solution to heat under reflux for 5 minutes. Crystals may begin to appear from the solution during this heating period.

Remove the heat source and allow the solution to cool slowly. The product will crystallize from the solution as it cools. When the solution has cooled to room temperature, collect the crystals by vacuum filtration using a small Büchner funnel. Rinse the crystals with three 2-mL portions of cold, glacial acetic acid. Allow the crystals to dry in the air overnight. Weigh the dried product and determine its melting point. Pure β -benzopinacolone melts at 182°C. Obtain the infrared spectrum using the dry film method (Technique 25, Section 25.4) or by an alternative method (Technique 25, Section 25.5) and the NMR spectrum in CDCl_3 (Technique 26, Section 26.1).

Calculate the percentage yield. Submit the product to your instructor in a labeled vial, along with your spectra. Interpret your spectra, showing how they are consistent with the rearranged structure of the product.

QUESTIONS

1. Can you think of a way to produce the benzophenone $n-\pi^*$ triplet T_1 *without* having benzophenone pass through its first singlet state? Explain.
2. A reaction similar to the one described here occurs when benzophenone is treated with the metal magnesium (pinacol reduction).

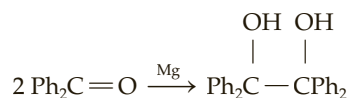

Compare the mechanism of this reaction with the photoreduction mechanism. What are the differences?

3. Which of the following molecules do you expect would be useful in quenching benzophenone photoreduction? Explain.

Oxygen	($S_1 = 22$ kcal/mol)
9,10-Diphenylanthracene	($T_1 = 42$ kcal/mol)
<i>trans</i> -1,3-Pentadiene	($T_1 = 59$ kcal/mol)
Naphthalene	($T_1 = 61$ kcal/mol)
Biphenyl	($T_1 = 66$ kcal/mol)
Toluene	($T_1 = 83$ kcal/mol)
Benzene	($T_1 = 84$ kcal/mol)

Fireflies and Photochemistry

The production of light as a result of a chemical reaction is called **chemiluminescence**. A chemiluminescent reaction generally produces one of the product molecules in an electronically excited state. The excited state emits a photon, and light is produced. If a reaction that produces light is biochemical, occurring in a living organism, the phenomenon is called **bioluminescence**.


The light produced by fireflies and other bioluminescent organisms has fascinated observers for many years. Many different organisms have developed the ability to emit light. They include bacteria, fungi, protozoans, hydras, marine worms, sponges, corals, jellyfish, crustaceans, clams, snails, squids, fish, and insects. Curiously, among the higher forms of life, only fish are included on the list. Amphibians, reptiles, birds, mammals, and the higher plants are excluded. Among the marine species, none is a freshwater organism. The excellent *Scientific American* article by McElroy and Seliger (see References) delineates the natural history, characteristics, and habits of many bioluminescent organisms.

The first significant studies of a bioluminescent organism were performed by the French physiologist Raphael Dubois in 1887. He studied the mollusk *Pholas dactylis*, a bioluminescent clam, indigenous to the Mediterranean Sea. Dubois found that a cold-water extract of the clam was able to emit light for several minutes following the extraction. When the light emission ceased, it could be restored, Dubois found, by a material extracted from the clam by hot water. A hot-water extract of the clam alone did not produce the luminescence. Reasoning carefully, Dubois concluded that there was an enzyme in the cold-water extract that was destroyed in hot water. The luminescent compound, however, could be extracted without destruction in either hot or cold water. He called the luminescent material **luciferin**, and the enzyme that induced it to emit light **luciferase**; both names were derived from *Lucifer*, a Latin name meaning “bearer of light.” Today the luminescent materials from all organisms are called *luciferins*, and the associated enzymes are called *luciferases*.


The most extensively studied bioluminescent organism is the firefly. Fireflies are found in many parts of the world and probably represent the most familiar example of bioluminescence. In such areas, on a typical summer evening, fireflies, or “lightning bugs,” can frequently be seen to emit flashes of light as they cavort over the lawn or in the garden. It is now universally accepted that the luminescence of fireflies is a mating strategy. The male firefly flies about 2 feet above the ground and emits flashes of light at regular intervals. The female, who remains stationary on the ground, waits a characteristic interval and then flashes a response. In return, the male reorients his direction of flight toward her and flashes a signal once again. The entire cycle is rarely repeated more than 5 to 10 times before the male reaches the female. Fireflies of different species can recognize one another by their flash patterns, which vary in number, rate, and duration among species.

Although the total structure of the luciferase enzyme of the American firefly *Photinus pyralis* is unknown, the structure of the luciferin has been established. In

spite of a large amount of experimental work, however, the complete nature of the chemical reactions that produce the light is still subject to some controversy. It is possible, nevertheless, to outline the most salient details of the reaction.

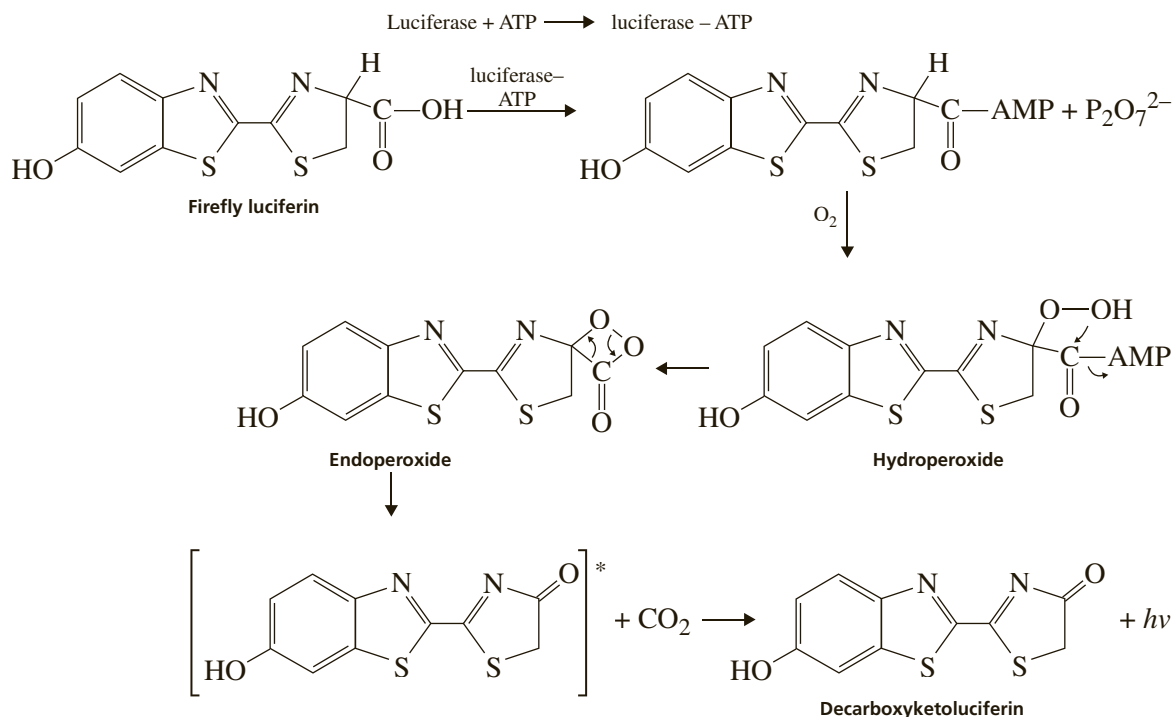

Besides luciferin and luciferase, other substances—magnesium(II), ATP (adenosine triphosphate), and molecular oxygen—are needed to produce the luminescence. In the postulated first step of the reaction, luciferase complexes with an ATP molecule. In the second step, luciferin binds to luciferase and reacts with the already-bound ATP molecule to become “primed.” In this reaction, pyrophosphate ion is expelled, and AMP (adenosine monophosphate) becomes attached to the carboxyl group of the luciferin. In the third step, the luciferin–AMP complex is oxidized by molecular oxygen to form a hydroperoxide; this cyclizes with the carboxyl group, expelling AMP and forming the cyclic endoperoxide. This reaction would be difficult if the carboxyl group of luciferin had not been primed with ATP. The endoperoxide is unstable and readily decarboxylates, producing decarboxylketoluciferin in an *electronically excited state*, which is deactivated by the emission of a photon (fluorescence). Thus, it is the cleavage of the four-membered-ring endoperoxide that leads to the electronically excited molecule and hence the bioluminescence.

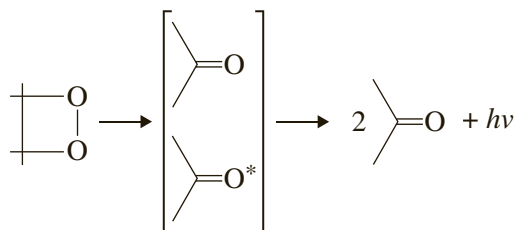


That one of the two carbonyl groups, either that of the decarboxylketoluciferin or that of the carbon dioxide, should be formed in an excited state can be readily predicted from the orbital symmetry conservation principles of Woodward and Hoffmann. This reaction is formally like the decomposition of a cyclobutane ring and yields two ethylene molecules. In analyzing the forward course of that

reaction, that is, 2 ethylene \longrightarrow cyclobutane, one can easily show that the reaction, which involves four π electrons, is forbidden for two ground-state ethylenes but allowed for only one ethylene in the ground state and the other in an excited state. This suggests that, in the reverse process, one of the ethylene molecules should be formed in an excited state. Extending these arguments to the endoperoxide also suggests that one of the two carbonyl groups should be formed in its excited state.

The emitting molecule, decarboxylketoluciferin, has been isolated and synthesized. When it is excited photochemically by photon absorption in basic solution ($\text{pH} > 7.5\text{--}8.0$), it fluoresces, giving a fluorescence emission spectrum that is identical to the emission spectrum produced by the interaction of firefly luciferin and firefly luciferase. The emitting form of decarboxylketoluciferin has thus been identified as the **enol dianion**. In neutral or acidic solution, the emission spectrum of decarboxylketoluciferin does not match the emission spectrum of the bioluminescent system.

The exact function of the enzyme firefly luciferase is not yet known, but it is clear that all these reactions occur while luciferin is bound to the enzyme as a substrate. Also, because the enzyme undoubtedly has several basic groups (—COO^- , —NH_2 , and so on), the buffering action of those groups would easily explain why the enol dianion is also the emitting form of decarboxylketoluciferin in the biological system.

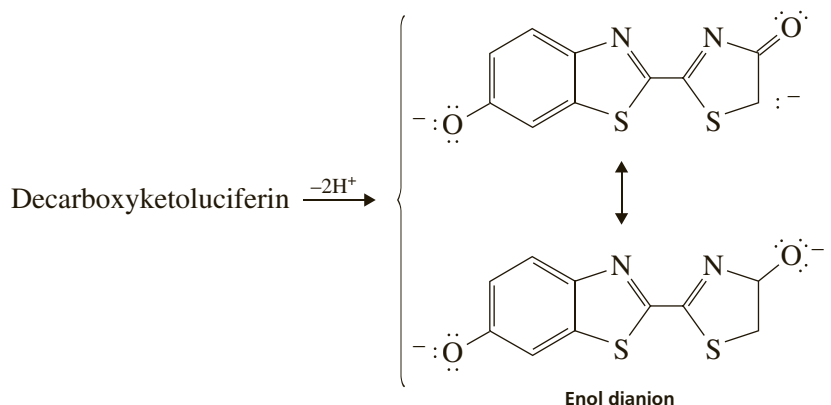

Most chemiluminescent and bioluminescent reactions require oxygen. Likewise, most produce an electronically excited emitting species through the decomposition of a **peroxide** of one sort or another. In the experiment that follows, a **chemiluminescent** reaction that involves the decomposition of a peroxide intermediate is described.

REFERENCES

- Clayton, R. K. The Luminescence of Fireflies and Other Living Things. In *Light and Living Matter. The Biological Part*; McGraw Hill: New York, 1971; Vol. 2.
- Fox, J. L. Theory May Explain Firefly Luminescence. *Chem. Eng. News*. **1978**, 56, 17.
- Harvey, E. N. *Bioluminescence*; Academic Press: New York, 1952.
- Hastings, J. W. Bioluminescence. *Annu. Rev. Biochem.* **1968**, 37, 597.
- McCapra, F. Chemical Mechanisms in Bioluminescence. *Acc. Chem. Res.* **1976**, 9, 201.
- McElroy, W. D.; Seliger, H. H. Biological Luminescence. *Sci. Am.* **1962**, 207, 76.
- McElroy, W. D.; Seliger, H. H.; White, E. H. Mechanism of Bioluminescence, Chemiluminescence and Enzyme Function in the Oxidation of Firefly Luciferin. *Photochem. Photobiol.* **1969**, 10, 153.
- Seliger, H. H.; McElroy, W. D. *Light: Physical and Biological Action*; Academic Press: New York, 1965.

Luminol


Chemiluminescence

Energy transfer

Reduction of a nitro group


Amide formation


In this experiment, the chemiluminescent compound **luminol**, or **5-amino-phthalhydrazide**, will be synthesized from 3-nitrophthalic acid.

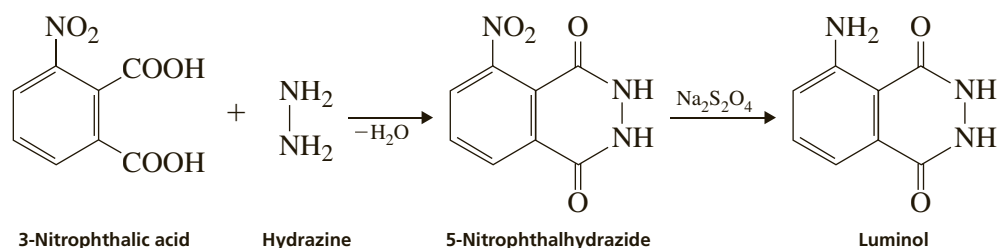


The first step of the synthesis is the simple formation of a cyclic diamide, 5-nitrophthalhydrazide, by reaction of 3-nitrophthalic acid with hydrazine. Reduction of the nitro group with sodium dithionite affords luminol.

In neutral solution, luminol exists largely as a dipolar anion (zwitterion). This dipolar ion exhibits a weak blue fluorescence after being exposed to light. However, in alkaline solution, luminol is converted to its dianion, which may be oxidized by molecular oxygen to give an intermediate that is chemiluminescent. The reaction is thought to have the following sequence:


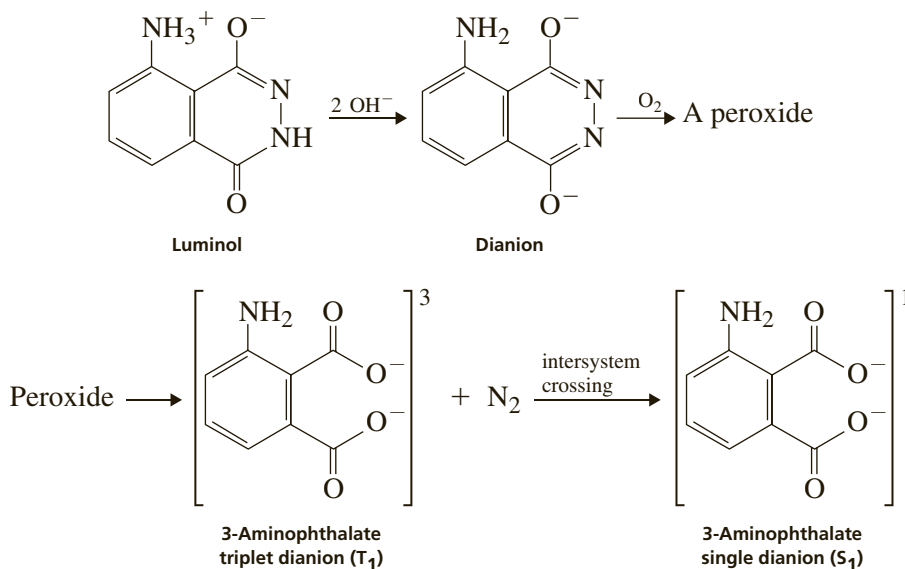


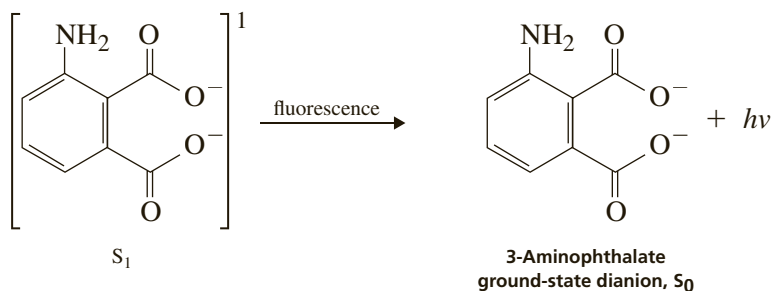

The dianion of luminol undergoes a reaction with molecular oxygen to form a peroxide of unknown structure. This peroxide is unstable and decomposes with the evolution of nitrogen gas, producing the 3-aminophthalate dianion in an electronically excited state. The excited dianion emits a photon that is visible as light. One very attractive hypothesis for the structure of the peroxide postulates a cyclic endoperoxide that decomposes by the following mechanism:

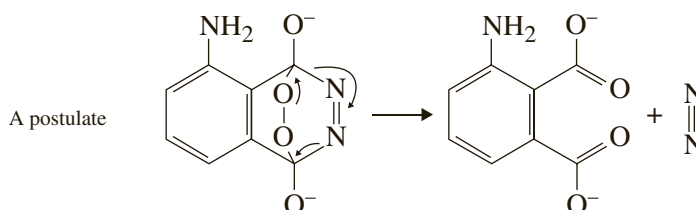

Certain experimental facts argue against this intermediate, however. For instance, certain acyclic hydrazides that cannot form a similar intermediate have also been found to be chemiluminescent.

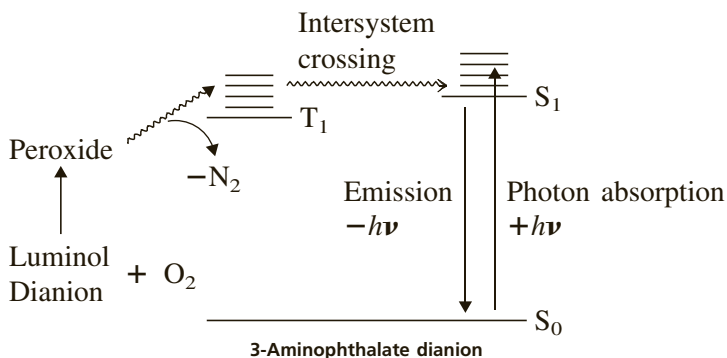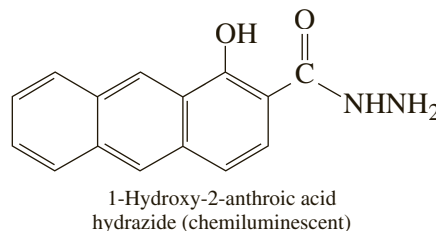

Fluorescence emission spectrum of the 3-aminophthalate dianion.

Although the nature of the peroxide is still debatable, the remainder of the reaction is well understood. The chemical products of the reaction have been shown to be the 3-aminophthalate dianion and molecular nitrogen. The intermediate that emits light has been identified definitely as the *excited state singlet* of the 3-aminophthalate dianion.¹ Thus, the fluorescence emission spectrum of the 3-aminophthalate dianion (produced by photon absorption) is identical to the spectrum of the light emitted from the chemiluminescent reaction. However, for numerous complicated reasons, it is believed that the

¹ The terms *singlet*, *triplet*, *intersystem crossing*, *energy transfer*, and *quenching* are explained in Experiment 55.

3-aminophthalate dianion is formed first as a vibrationally excited triplet state molecule, which makes the intersystem crossing to the singlet state before emission of a photon.

The excited state of the 3-aminophthalate dianion may be quenched by suitable acceptor molecules, or the energy (about 50–80 Kcal/mol) may be transferred to give emission from the acceptor molecules. Several such experiments are described in the following procedure.

The system chosen for the chemiluminescence studies of luminol in this experiment uses dimethylsulfoxide (CH_3)₂SO as the solvent, potassium hydroxide as the base required for the formation of the dianion of luminol, and molecular oxygen. Several alternative systems have been used, substituting hydrogen peroxide and an oxidizing agent for molecular oxygen. An aqueous system using potassium ferricyanide and hydrogen peroxide is an alternative system used frequently.

REFERENCES

- Rahaut, M. M. Chemiluminescence from Concerted Peroxide Decomposition Reactions. *Acc. Chem. Res.*, **1969**, *2*, 80.
- White, E. H.; Roswell, D. F. The Chemiluminescence of Organic Hydrazides. *Acc. Chem. Res.*, **1970**, *3*, 54.

REQUIRED READING

- Review:* Technique 3 Reaction Methods, Section 7.10
- New:* Essay Fireflies and Photochemistry

SPECIAL INSTRUCTIONS

This entire experiment can be completed in about 1 hour. When you are working with hydrazine, you should remember that it is toxic and should not be spilled on the skin. It is also a suspected carcinogen. Dimethylsulfoxide may also be toxic; avoid breathing the vapors or spilling it on your skin.

A darkened room is required to observe adequately the chemiluminescence of luminol. A darkened hood that has had its window covered with butcher paper or aluminum foil also works well. Other fluorescent dyes besides those mentioned (for instance, 9,10-diphenylanthracene) can also be used for the energy-transfer experiments. The dyes selected may depend on what is immediately available. The instructor may have each student use one dye for the energy-transfer experiments, with one student making a comparison experiment without a dye.

SUGGESTED WASTE DISPOSAL

Dispose of the filtrate from the vacuum filtration of 5-nitrophthalhydrazide in the container designated for nonhalogenated organic solvents. The filtrate from the vacuum filtration of 5-aminophthalhydrazide may be diluted with water and

poured into the waste container designated for aqueous waste. The mixture containing potassium hydroxide, dimethylsulfoxide, and luminol should be placed in the special container designated for this material.

PROCEDURE

Part A. 5-Nitrophthalhydrazide

Place 0.300 g of 3-nitrophthalic acid and 0.4 mL of a 10% aqueous solution of hydrazine (use gloves) in a small side-arm test tube.² At the same time, heat 4 mL of water in a beaker on a hot plate to about 80°C. Heat the test tube over a microburner until the solid dissolves. Add 0.8 mL of triethylene glycol and clamp the test tube in an upright position on a ring stand. Place a thermometer (do not seal the system) and a boiling stone in the test tube and attach a piece of pressure tubing to the side arm. Connect this tubing to an aspirator (use a trap). The thermometer bulb should be in the liquid as much as possible. Heat the solution with a microburner until the liquid boils vigorously and the refluxing water vapor is drawn away by the aspirator vacuum (the temperature will rise to about 120°C). Continue heating and allow the temperature to increase rapidly until it rises just above 200°C. This heating requires 1–2 minutes, and you must watch the temperature closely to avoid heating the mixture well above 200°C. Remove the burner briefly when this temperature has been achieved and then resume gentle heating to maintain a fairly constant temperature of 210–220°C for about 2 minutes. Allow the test tube to cool to about 100°C, add the 4.0 mL of hot water that was prepared previously, and cool the test tube to room temperature by allowing tap water to flow over the outside of the test tube. Collect the brown crystals of 5-nitrophthalhydrazide by vacuum filtration, using a small Hirsch funnel. It is not necessary to dry the product before you go on with the next reaction step.

Part B. Luminol (5-Aminophthalhydrazide)

Transfer the moist 5-nitrophthalhydrazide to a 13 x 100-mm test tube. Add 1.30 mL of a 10% sodium hydroxide solution and agitate the mixture until the hydrazide dissolves. Add 0.80 g of sodium dithionite dihydrate (sodium hydrosulfite dihydrate, $\text{Na}_2\text{S}_2\text{O}_4 \cdot 2\text{H}_2\text{O}$). Using a Pasteur pipette, add 1–2 mL of water to wash the solid from the walls of the test tube. Add a boiling stone to the test tube. Heat the test tube until the solution boils, agitate the solution, and maintain the boiling, continuing agitation, for 5 minutes. Add 0.50 mL of glacial acetic acid and cool the test tube to room temperature by allowing tap water to flow over the outside of it. Agitate the mixture during the cooling step. Collect the light yellow or gold crystals of luminol by vacuum filtration, using a small Hirsch funnel. Save a small sample of this product, allow it to dry overnight, and determine its melting point (mp 319–320°C). The remainder of the luminol may be used without drying for the chemiluminescence experiments.

CAUTION



Be careful not to let any of the mixture touch your skin while shaking the flask. Hold the stopper securely.

Part C. Chemiluminescence Experiments

Cover the bottom of a 10-mL Erlenmeyer flask with a layer of potassium hydroxide pellets. Add enough dimethylsulfoxide to cover the pellets. Add about 0.025 g of the moist luminol to the flask, stopper it, and shake it vigorously to mix

² A 10% aqueous solution of hydrazine can be prepared by diluting 15.6 g of a commercial 64% hydrazine solution to a volume of 100 mL using water.

air into the solution.³ In a dark room, a faint glow of bluish white light will be visible. The intensity of the glow will increase with continued shaking of the flask and occasional removal of the stopper to admit more air.

To observe energy transfer to a fluorescent dye, dissolve one or two crystals of the indicator dye in about 0.25 mL of water. Add the dye solution to the dimethylsulfoxide solution of luminol, stopper the flask, and shake the mixture vigorously. Observe the intensity and the color of the light produced.

The following table shows some dyes and the colors produced when they are mixed with luminol. Other dyes not included on this list may also be tested in this experiment.

<i>Fluorescent Dye</i>	<i>Color</i>
No dye	Faint blush white
2,6-Dichloroindophenol	Blue
9-Aminoacridine	Blue-green
Eosin	Salmon pink
Fluorescein	Yellow-green
Dichlorofluorescein	Yellow-orange
Rhodamine B	Green
Phenolphthalein	Purple

³An alternative method for demonstrating chemiluminescence, using potassium ferricyanide and hydrogen peroxide as oxidizing agents, is described in E. H. Huntress, L. N. Stanley, and A. S. Parker, *Journal of Chemical Education*, 11 (1934): 142.

Identification of Organic Substances

Identification of Unknowns


Qualitative organic analysis, the identification and characterization of unknown compounds, is an important part of organic chemistry. Every chemist must learn the appropriate methods for establishing the identity of a compound. In this experiment, you will be issued an unknown compound and will be asked to identify it through chemical and spectroscopic methods. Your instructor may give you a general unknown or a specific unknown. With a **general unknown**, you must first determine the class of compound to which the unknown belongs, that is, identify its main functional group; then you must determine the specific compound in that class that corresponds to the unknown. With a **specific unknown**, you will know the class of compound (ketone, alcohol, amine, and so on) in advance, and it will be necessary to determine only whatever specific member of that class was issued to you as an unknown. This experiment is designed so that the instructor can issue several general unknowns or as many as six successive specific unknowns, each having a different main functional group.

Although there are millions of organic compounds that an organic chemist might be called on to identify, the scope of this experiment is necessarily limited. In this textbook, about 500 compounds are included in the tables of possible unknowns given for the experiment (see Appendix 1). Your instructor may wish to expand the list of possible unknowns, however. In such a case, you will have to consult more extensive tables, such as those found in the work compiled by Rapoport (see References). In addition, the experiment is restricted to include only seven important functional groups:


Aldehydes	Amines
Ketones	Alcohols
Carboxylic acids	Esters
Phenols	

Even though this list of functional groups omits some of the important types of compounds (alkyl halides, alkenes, alkynes, aromatics, ethers, amides, mercaptans, nitriles, acid chlorides, acid anhydrides, nitro compounds, and so on), the methods introduced here can be applied equally well to other classes of compounds. The list is sufficiently broad to illustrate all the principles involved in identifying an unknown compound.

In addition, although many of the functional groups listed as being excluded will not appear as the major functional group in a compound, several of them will frequently appear as secondary, or subsidiary, functional groups. Three examples of this are presented here.

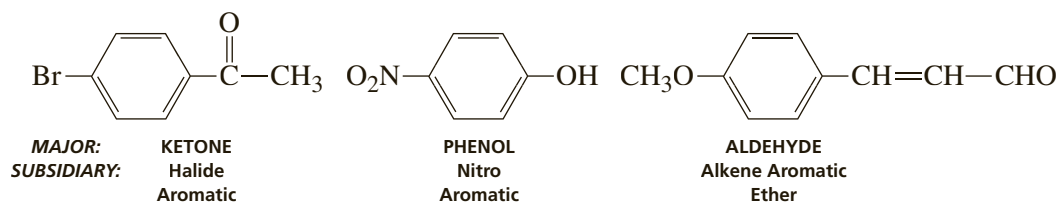

The groups included that have subsidiary status are

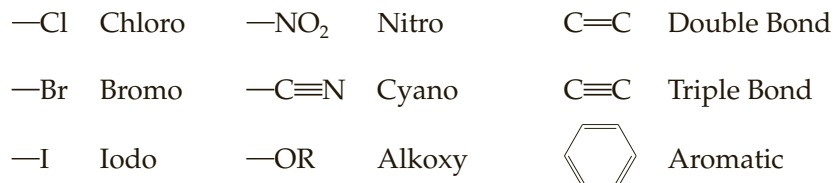

The experiment presents all of the chief chemical and spectroscopic methods of determining the main functional groups, and it includes methods for verifying the presence of the subsidiary functional groups as well. It will usually not be necessary to determine the presence of the subsidiary functional groups to identify the unknown compound correctly. Every piece of information helps the identification, however, and if these groups can be detected easily, you should not hesitate to determine them. Finally, complex bifunctional compounds are generally avoided in this experiment; only a few are included.

PROCEDURE

How to Proceed— Option 1

Fortunately, we can detail a fairly straightforward procedure for determining all of the necessary pieces of information. This procedure consists of the following steps:

Part One: Chemical Classification

1. Preliminary classification by physical state, color, and odor
2. Melting-point or boiling-point determination; other physical data
3. Purification, if necessary
4. Determination of solubility behavior in water and in acids and bases
5. Simple preliminary tests: Beilstein, ignition (combustion)
6. Application of relevant chemical classification tests
7. Inspection of tables for possible structure(s) of unknown; elimination of unlikely compounds

Part Two: Spectroscopy

8. Determination of infrared and NMR spectra

Part Three: Optional Procedures

9. Elemental analysis, if necessary
10. Preparation of derivatives, if required
11. Confirmation of identity

Each of these steps is discussed briefly starting below:

**Green Chemistry
Method: How to
Proceed—Option 2**

At the option of your instructor, another approach may be taken in determining the structure of unknowns in the organic laboratory. This approach makes minimal use of classification tests but retains the solubility tests as the main way of determining functional groups and spectroscopy as a way of determining the detailed structure of an unknown. Elimination of classification tests described in Part One, number 6, tremendously reduces the waste generated in the laboratory. It also eliminates the use of many of the toxic and potentially dangerous reagents that are a standard part of the traditional classification tests. This approach is, therefore, a “Green” approach to solving structures of organic compounds.

Although classification tests can be useful in determining the identity of an unknown compound, spectroscopic methods have become the principal means by which an organic chemist identifies unknown substances. The technology and instrumentation available has almost obviated the need for classification tests, because valuable information can be discovered simply by obtaining infrared and NMR spectra. Option 2 relies heavily on the spectroscopic results; if acetone- d_6 or DMSO- d_6 are used as NMR spectroscopy solvents, this becomes a more environmentally sound approach.

The ability to use IR and NMR spectroscopy and evaluate spectra inherently requires a logical sequence of steps in the identification of an unknown. By relying on these techniques, students learn the techniques and higher-order thinking skills that they would be required to know and use for a career in chemistry. This approach more closely simulates the types of structure-proof methods that one would find in a modern research or industrial laboratory. Students can still learn how to go through the logical steps used in the classification tests by practicing these methods in a more environmentally friendly scenario through the use of computer simulations.

The procedure for determining the structure of a compound using the environmentally friendly approach is fairly straightforward and consists of the following steps:

Part One: Chemical Classification

1. Preliminary classification by physical state, color, and odor
2. Melting-point or boiling-point determination; other physical data
3. Purification, if necessary
4. Determination of solubility behavior in water and in acids and bases
5. Simple preliminary tests: Beilstein, ignition (combustion)
6. Inspection of tables for possible structure(s) of unknowns

Part Two: Spectroscopy

7. Determination of infrared and NMR (proton and ^{13}C , if available) spectra
8. Confirmation of structure

In many cases, the type of compound and functional group should be discovered after completing Part One. Spectroscopy (Part Two) will be used *principally* to confirm the structural assignment and to provide further information toward identifying the unknown. Your instructor may not allow you to obtain spectroscopic information (infrared or NMR) until you have completed Part One. Show your test results to your instructor for approval. Once this part has been completed, you should have narrowed the list of possible compounds to a few likely candidates, *all containing the same functional group*. In other words, you should have determined the principal functional group. You *must* obtain approval from the instructor to perform spectroscopy.

The functional groups that may be included in the unknowns are listed on the first page of this experiment. Tables of possible compounds are listed in the Appendix 1 of this book.

1. PRELIMINARY CLASSIFICATION

Note the physical characteristics of the unknown, including its color, odor, and physical state (liquid, solid, crystalline form). Many compounds have characteristic colors or odors, or they crystallize with a specific crystal structure. This information can often be found in a handbook and can be checked later. Compounds with a high degree of conjugation are frequently yellow to red. Amines often have a fishlike odor. Esters have a pleasant fruity or floral odor. Acids have a sharp and pungent odor. A part of the training of every good chemist includes cultivating the ability to recognize familiar or typical odors. As a note of caution, many compounds have distinctly unpleasant or nauseating odors. Some have corrosive vapors. Sniff any unknown substance with the greatest caution. As a first step, open the container, hold it away from you, and using your hand, carefully waft the vapors toward your nose. If you get past this stage, a closer inspection will be possible.

2. MELTING-POINT OR BOILING-POINT DETERMINATION

The single most useful piece of information to have for an unknown compound is its melting point or boiling point. Either piece of data will drastically limit the compounds that are possible. The electric melting-point apparatus gives a rapid and accurate measurement (see Technique 9, Sections 9.5 and 9.7). To save time, you can often determine two separate melting points. The first determination can be made rapidly to get an approximate value. Then you can determine the second melting point more carefully. Because some of the unknown solids contain traces of impurities, you may find that your observed melting point is lower than the values found in the tables in Appendix 1. This is especially true for low-melting compounds ($<50^{\circ}\text{C}$). For these low-melting compounds, it is a good idea to look at compounds in the tables in Appendix 1 that have melting points above your observed melting-point range. The same advice may apply to other solid compounds issued to you as unknowns.

The boiling point is easily obtained by a simple distillation of the unknown (see Technique 14, Section 14.3) by reflux (see Technique 13, Section 13.2), by a microboiling-point determination (see Technique 13, Section 13.2), or by Vernier LabPro interface method (see Technique 13, Section 13.5). The simple distillation has the advantage in that it also purifies the compound. The smallest distilling flask available should be used if a simple distillation is performed, and you should be sure that the thermometer bulb is fully immersed in the vapor of the distilling liquid. The liquid should be distilled rapidly to determine an accurate boiling-point value. The microboiling-point method requires the least amount of unknown, but the refluxing method is more reliable and requires much less liquid than that required for distillation.

When inspecting the tables of unknowns in Appendix 1, you may find that the observed boiling point that you determined is lower than the value for the corresponding compound listed in the tables. This is especially true for compounds boiling above 200°C . It is less likely, but not impossible, that the observed boiling point of your unknown will be higher than the value given in the table. Thus, your strategy should be to look for boiling points of compounds in the tables that are

nearly equal to or above the value you obtained, within a range of about $\pm 5^\circ\text{C}$. For high-boiling liquid compounds ($>200^\circ\text{C}$), you may need to apply a thermometer correction (see Technique 13, Section 13.3).

3. PURIFICATION

If the melting point of a solid has a wide range (about 5°C), the solid should be recrystallized and the melting point redetermined.

If a liquid was highly colored before distillation, if it yielded a wide boiling-point range, or if the temperature did not hold constant during the distillation, it should be redistilled to determine a new temperature range. A reduced-pressure distillation is in order for high-boiling liquids or for those that show any sign of decomposition on heating.

Occasionally, column chromatography may be necessary to purify solids that have large amounts of impurities and do not yield satisfactory results on crystallization.

Acidic or basic impurities that contaminate a neutral compound may often be removed by dissolving the compound in a low-boiling solvent, such as CH_2Cl_2 or ether, and extracting with 5% NaHCO_3 or 5% HCl , respectively. Conversely, acidic or basic compounds can be purified by dissolving them in 5% NaHCO_3 or 5% HCl , respectively, and extracting them with a low-boiling organic solvent to remove impurities. After the aqueous solution has been neutralized, the desired compound can be recovered by extraction.

4. SOLUBILITY BEHAVIOR

Tests on solubility are described fully in Experiment 57A. They are extremely important. Determine the solubility of small amounts of the unknown in water, 5% HCl , 5% NaHCO_3 , 5% NaOH , concentrated H_2SO_4 , and organic solvents. This information reveals whether a compound is an acid, a base, or a neutral substance. The sulfuric acid test reveals whether a neutral compound has a functional group that contains an oxygen, a nitrogen, or a sulfur atom that can be protonated. This information allows you to eliminate or to choose various functional-group possibilities. The solubility tests must be made on *all* unknowns. It may be helpful to consult the Merck Index for the solubility of your compound in organic solvents. Checking the solubility of the compound in these solvents can sometimes help to confirm the identity of your compound.

5. PRELIMINARY TESTS

The two combustion tests, the Beilstein test (Experiment 57B) and the ignition test (Experiment 57C), can be performed easily and quickly, and they often give valuable information. It is recommended that they be performed on all unknowns.

6. CHEMICAL CLASSIFICATION TESTS

The solubility tests usually suggest or eliminate several possible functional groups. The chemical classification tests listed in Experiment 57 allow you to distinguish among the possible choices. Choose only those tests that the solubility

tests suggest might be meaningful. Time will be wasted performing unnecessary tests. There is no substitute for a firsthand, thorough knowledge of these tests. Study each of the sections carefully until you understand the significance of each test. Also, it is essential to actually try the tests on *known* substances. In this way, it will be easier to recognize a positive test. Appropriate test compounds are listed for many of the tests. When you are performing a test that is new to you, it is always good practice to run the test separately on both a known substance and the unknown *at the same time*. This practice lets you compare results directly.

Do not perform the chemical tests either haphazardly or in a methodical, comprehensive sequence. Instead, use the tests selectively. Solubility tests automatically eliminate the need for some of the chemical tests. Each successive test will either eliminate the need for another test or dictate its use. You should also examine the Tables of Unknowns in Appendix 1 carefully. The boiling point or the melting point of the unknown may eliminate the need for many of the tests. For instance, the possible compounds may simply not include one with a double bond. *Efficiency* is the key word here. Do not waste time performing nonsensical or unnecessary tests. Many possibilities can be eliminated on the basis of logic alone.

How you proceed with the following steps may be limited by your instructor's wishes. Many instructors may restrict your access to infrared and NMR spectra until you have narrowed your choices to a few compounds, *all within the same class*. Others may have you determine these data routinely. Some instructors may want students to perform elemental analysis on all unknowns; others may restrict it to only the most essential situations. Again, some instructors may require derivatives as a final confirmation of the compound's identity; others may not wish to use them at all.

7. INSPECTION OF TABLES FOR POSSIBLE STRUCTURES

Once the melting or boiling point, the solubilities, and the main chemical classification tests have been made, you should be able to identify the class of compound (aldehyde, ketone, and so on). At this stage, with the melting point or boiling point as a guide, you can compile a list of possible compounds from one of the appropriate tables in Appendix 1. It is very important to draw out the structures of compounds that fit the solubility, classification tests, and melting point or boiling point that were determined. If necessary, you can look up the structures in the *CRC Handbook*, *The Merck Index*, or the *Aldrich Handbook*. Remember that the boiling point or melting point recorded in the table may be higher than what you obtained in the laboratory (see Section 2 above).

The short list that you developed by inspection of the tables in Appendix 1 and the structures drawn should suggest that some additional tests may be needed to distinguish among the possibilities. For instance, one compound may be a methyl ketone, and the other may not. The iodoform test is called for to distinguish the two possibilities. The tests for the subsidiary functional groups may also be required. These tests are described in Experiments 57B and 57C. These tests should also be studied carefully; there is no substitute for firsthand knowledge about these tests.

8. SPECTROSCOPY

Spectroscopy is probably the most powerful and modern tool available to the chemist for determining the structure of an unknown compound. It is often possible to determine the structure through spectroscopy alone. On the other hand, there are

also situations for which spectroscopy may not be of much help, and the traditional methods must be relied on. For this reason, you should not use spectroscopy to the exclusion of the more traditional tests but rather as a confirmation of those results. Nevertheless, the main functional groups and their immediate environmental features can be determined quickly and accurately with spectroscopy.

9. ELEMENTAL ANALYSIS

Elemental analysis—which allows you to determine the presence of nitrogen, sulfur, or a specific halogen atom (Cl, Br, I) in a compound—is often useful; however, other information may render these tests unnecessary. A compound identified as an amine by solubility tests obviously contains nitrogen. Many nitrogen-containing groups (for instance, nitro groups) can be identified by infrared spectroscopy. Finally, it is not usually necessary to identify a specific halogen. The simple information that the compound contains a halogen (any halogen) may be enough information to distinguish between two compounds. A simple Beilstein test provides this information.

10. DERIVATIVES

One of the principal tests for the correct identification of an unknown compound is to convert the compound by a chemical reaction to another known compound. This second compound is called a **derivative**. The best derivatives are solid compounds, because the melting point of a solid provides an accurate and reliable identification of most compounds. Solids are also easily purified through crystallization. The derivative provides a way of distinguishing two otherwise very similar compounds. Usually, they will have derivatives (both prepared by the same reaction) that have different melting points. Tables of unknowns and derivatives are listed in Appendix 1. Procedures for preparing derivatives are given in Appendix 2.

11. CONFIRMATION OF IDENTITY

A rigid and final test for identifying an unknown can be made if an “authentic” sample of the compound is available for comparison. One can compare infrared and NMR spectra of the unknown compound with the spectra of the known compound. If the spectra match, peak for peak, then the identity is probably certain. Other physical and chemical properties can also be compared. If the compound is a solid, a convenient test is the mixture melting point (see Technique 9, Section 9.4). Thin-layer or gas-chromatographic comparisons may also be useful. For thin-layer analysis, however, it may be necessary to experiment with several different development solvents to reach a satisfactory conclusion about the identity of the substance in question.

Although we cannot be complete in this experiment in terms of the functional groups covered or the tests described, the experiment should provide a good introduction to the methods and the techniques chemists use to identify unknown compounds. Textbooks that cover the subject more thoroughly are listed in the References. You are encouraged to consult these for more information, including specific methods and classification tests.

REFERENCES

Comprehensive Textbooks

- Cheronis, N. D.; Entrikin, J. B. *Identification of Organic Compounds*; Wiley-Interscience: New York, 1963.
- Pasto, D. J.; Johnson, C. R. *Laboratory Text for Organic Chemistry*; Prentice Hall: Englewood Cliffs, NJ, 1979.
- Shriner, R. L.; Hermann, C. K. F.; Morrill, T. C.; Curtin, D. Y.; Fuson, R. C. *The Systematic Identification of Organic Compounds*, 8th ed.; Wiley: New York, 2004.

Spectroscopy

- Bellamy, L. J. *The Infra-red Spectra of Complex Molecules*, 3rd ed.; Methuen: New York, 1975.
- Colthup, N. B.; Daly, L. H.; Wiberly, S. E. *Introduction to Infrared and Raman Spectroscopy*, 3rd ed.; Academic Press: San Diego, CA, 1990.
- Lin-Vien, D.; Colthup, N. B.; Fateley, W. B.; Grasselli, J. G. *The Handbook of Infrared and Raman Characteristic Frequencies of Organic Molecules*; Academic Press: San Diego, CA, 1991.
- Nakanishi, K.; Solomon, P.H. *Infrared Absorption Spectroscopy*, 2nd ed.; Holden-Day: San Francisco, 1998.
- Pavia, D. L.; Lampman, G. M.; Kriz, G. S.; Vyryan, J. R. *Introduction to Spectroscopy: A Guide for Students of Organic Chemistry*, 5th ed.; Cengage Learning: Stamford, CT, 2015.
- Silverstein, R. M.; Webster, F. X.; Kiemle, D. *Spectrometric Identification of Organic Compounds*, 7th ed.; Wiley: New York, 2005.

Extensive Tables of Compounds and Derivatives

- Rappoport, Z., Ed. *Handbook of Tables for Organic Compound Identification*, 3rd ed.; CRC Press: Boca Raton, FL, 1967.


57A EXPERIMENT 57A

Solubility Tests

Solubility tests should be performed on *every unknown*. They are extremely important in determining the nature of the main functional group of the unknown compound. Although the tests are very simple, it can take a little practice to interpret the results correctly. Therefore, it is important that you follow the instructions carefully. It is also very important that you determine the solubility of the five test compounds that will be available as unknowns. You should verify your results with your instructor. Solubility tests reveal whether the compound is a strong base (amine), a weak acid (phenol), a strong acid (carboxylic acid), or a neutral substance (aldehyde, ketone, alcohol, ester). The common solvents used to determine solubility types are

5% HCl	Concentrated H ₂ SO ₄
5% NaHCO ₃	Water
5% NaOH	Organic solvents

The solubility chart given in the next page indicates solvents in which compounds containing the various functional groups are likely to dissolve. The summary charts in Experiments 57D through 57I repeat this information for each functional group included in this experiment. In this section, the correct procedure for determining whether a compound is soluble in a test solvent is given. Also given is a series of explanations detailing the reasons that compounds having specific functional groups are soluble only in specific solvents. This is accomplished by indicating the type of chemistry or the type of chemical interaction that is possible in each solvent.

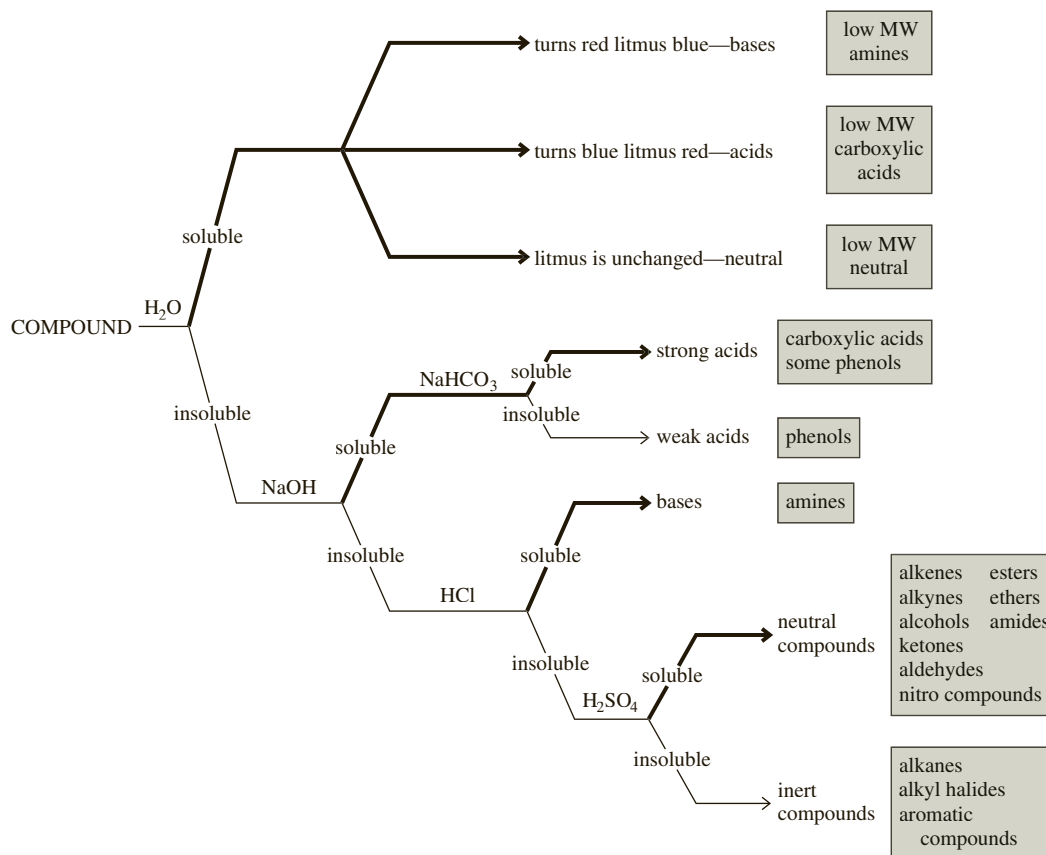

Solubility chart for compounds containing various functional groups. The most common functional groups in each class are printed in **bold-face type**.

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous solutions in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container.

SOLUBILITY TESTS

Procedure

Place about 2 mL of the solvent in a small test tube. Add 3 drops of an unknown liquid from a Pasteur pipette or a small amount (estimated at about 20 mg or 0.030 g) of an unknown solid into the solvent. Vigorously agitate the contents of the test tube by striking the side of the test tube rapidly with your fingers. Alternatively, use the rounded end of a microspatula (the larger style shown in Figure 10 in Experiment 1, "Introduction to Microscale Laboratory") to stir each sample by twirling the spatula rapidly. With either method, you should continue to mix the sample for about 3 minutes. After this time, note whether the compound is *soluble* (dissolves completely) or *insoluble* (none or very little of it dissolves). If the sample dissolves quickly, it may be helpful to determine the extent of solubility by adding several more drops of the liquid or another 20 mg of the solid. You should also stir this mixture for about 3 minutes. Solids in the form of large crystals need more time to dissolve than powders or very small crystals.

In some cases, it is helpful to crush the large crystals. Sometimes gentle heating helps, but strong heating is discouraged. When colored compounds dissolve, the solution assumes the color. When testing liquid unknowns, you should watch closely the upper portion of the solvent as you add the unknown. If mixing lines appear in the solvent, this is a good indication that the unknown is soluble. If the unknown liquid is insoluble, you may be able to see droplets of the unknown throughout the solvent when you agitate the mixture.

Using the procedure above and the Solubility chart, determine the solubility class of all unknown compounds issued to you. The solvents that will be most helpful in determining the solubility class are water, 5% HCl, 5% NaOH, and 5% NaHCO₃. You should always start at the left side of the Solubility chart and then move to the right, depending on whether or not your unknown is soluble. For example, if your unknown is insoluble in water, you would then test the solubility in NaOH. Note that most organic compounds will be only slightly soluble or insoluble in water. You need not use all of the solvents shown in this chart with every unknown. For example, if you find that the unknown is soluble in NaOH, then try NaHCO₃. If you find that it is soluble in *both* NaOH and NaHCO₃ then the unknown is likely to be a carboxylic acid and you can stop at this point. If the unknown is soluble in NaOH but not NaHCO₃ then the unknown is likely to be a phenol. Concentrated H₂SO₄ is always the last choice since it will “dissolve” or change the appearance of most organic compounds, except inert compounds. Sulfuric acid may yield a color change rather than dissolving the compound. A color change should be regarded as a positive solubility test. For the unknowns suggested in this textbook, checking for solubility in sulfuric acid will usually not be required. You should also keep in mind that some unknowns may have borderline solubility characteristics and may not behave exactly as indicated in the Solubility chart.

The typical solubility classes and conclusions are as follows:

Strong acids¹: soluble in NaHCO₃ and NaOH (carboxylic acids)

Weak Acids: soluble in NaOH, insoluble in NaHCO₃ (phenols)

Bases: soluble in HCl (amines)

Neutral compounds suggested in this book as unknowns (alcohols, ketones, aldehydes, and esters)

Organic compounds will generally be soluble in an organic solvent. The exceptions may be sodium or potassium salts of carboxylic acids or amine salts. These compounds will be soluble in water, and insoluble in solvents such as ether or dichloromethane (methylene chloride).

If a compound is found to dissolve in water, the pH of the aqueous solution should be estimated with pH paper or litmus. Compounds soluble in water are usually soluble in *all* the aqueous solvents. If a compound is only slightly soluble in water, it may be *more* soluble in another aqueous solvent. For instance, carboxylic acid may be only slightly soluble in water but very much soluble in dilute base.

Test Compounds

Five solubility unknowns will be used as practice to help you make proper observations and conclusions. This set of compounds can be found on the supply

¹ In General Chemistry, the term “strong acid” refers to acids such as HCl and H₂SO₄ that are completely ionized in water. In this book, we use the term “strong acid” to identify carboxylic acids and some phenols. Phenols that are very weak are identified as being “weak acids.”

shelf. The five unknowns include a strong acid (carboxylic acid), a weak acid (phenol), a base (amine), a neutral substance with an oxygen-containing functional group, and a substance that is inert. The general approach is to start at the left of the Solubility chart shown near the beginning of this section. Most organic compounds are either sparingly soluble in water or insoluble in water. If you find that the unknown is insoluble or only slightly soluble in water, then you should try aqueous sodium hydroxide, proceeding toward the right in the chart until you reach a conclusion. Sometimes you may find that a compound is slightly soluble in one reagent, but much more soluble in another. In that case, you should conclude that the unknown resides in the solubility class where it is *most soluble*. Solubility in sulfuric acid should always be the last solvent selected. Assign each of the unknowns in solubility classes, as follows:

- Carboxylic acid: soluble in NaOH and NaHCO₃
- Phenol: soluble in NaOH, but insoluble in NaHCO₃
- Amine: soluble in HCl, but insoluble in both NaOH and NaHCO₃
- Neutral compounds, ketone, aldehyde, alcohol, esters: insoluble in all except sulfuric acid
- Inert compounds (see chart for examples): insoluble in all, including sulfuric acid

This set of five unknowns will help you make proper observations. They will not be used further in your organic laboratory course. Note that almost all organic compounds, except inert ones, will be soluble in sulfuric acid. This reagent should always be the last one tried, as shown in the Solubility chart. Using solubility tests, distinguish these unknowns by type. Verify your answer with the instructor. Read the discussion sections that follow for details on solubility behavior. A more general discussion of solubility behavior is provided in Technique 10, Section 10.2.


Discussion

Solubility in Water

Compounds that contain four or fewer carbons and also contain oxygen, nitrogen, or sulfur are often soluble in water. Almost any functional group containing these elements will lead to water solubility for low-molecular-weight (C₄) compounds. Compounds having five or six carbons and any of those elements are often insoluble in water or have borderline solubility. Branching of the alkyl chain in a compound lowers the intermolecular forces between its molecules. This is usually reflected in a lowered boiling point or melting point and a greater solubility in water for the branched compound than for the corresponding straight-chain compound. This occurs simply because the molecules of the branched compound are more easily separated from one another. Thus, *t*-butyl alcohol would be expected to be more soluble in water than *n*-butyl alcohol.


When the ratio of the oxygen, nitrogen, or sulfur atoms in a compound to the carbon atoms is increased, the solubility of that compound in water often increases. This is due to the increased number of polar functional groups. Thus, 1,5-pentanediol would be expected to be more soluble in water than 1-pentanol.

As the size of the alkyl chain of a compound is increased beyond about four carbons, the influence of a polar functional group is diminished, and the water solubility begins to decrease. A few examples of these generalizations are given here.

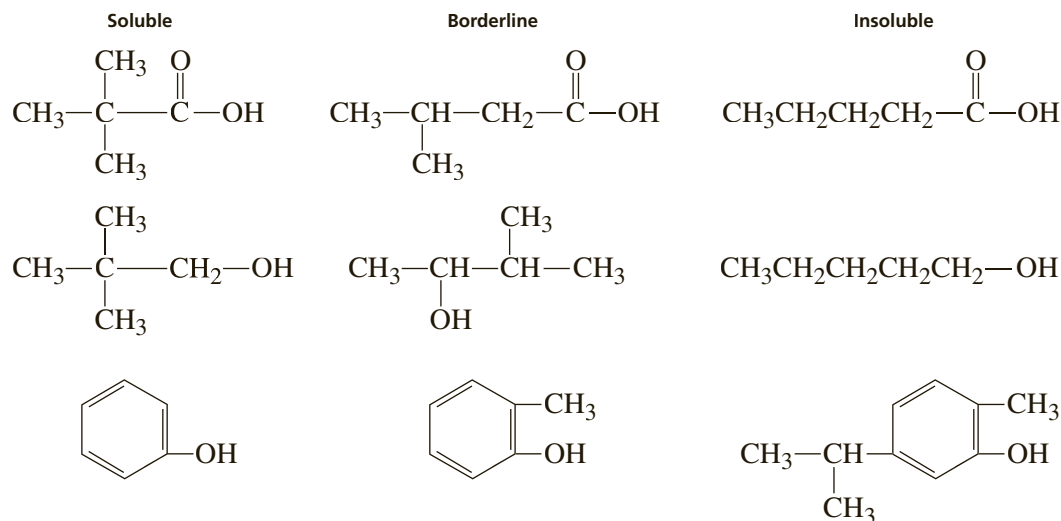


Solubility in 5% HCl

The possibility of an amine should be considered immediately if a compound is soluble in dilute acid (5% HCl). Aliphatic amines (RNH_2 , R_2NH , R_3N) are basic compounds that readily dissolve in acid because they form hydrochloride salts that are soluble in the aqueous medium:

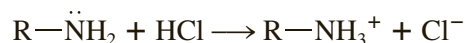

The substitution of an aromatic (benzene) ring Ar for an alkyl group R reduces the basicity of an amine somewhat, but the amine will still protonate, and it will still generally be soluble in dilute acid. The reduction in basicity in an aromatic amine is due to the resonance delocalization of the unshared electrons on the amino nitrogen of the free base. The delocalization is lost on protonation, a problem that does not exist for aliphatic amines. The substitution of two or three aromatic rings on an amine nitrogen reduces the basicity of the amine even further. Diaryl and triaryl amines do not dissolve in dilute HCl because they do not protonate easily. Thus, Ar_2NH and Ar_3N are insoluble in dilute acid. Some amines of very high molecular weight, such as tribromoaniline ($\text{MW} = 330$), may also be insoluble in dilute acid.

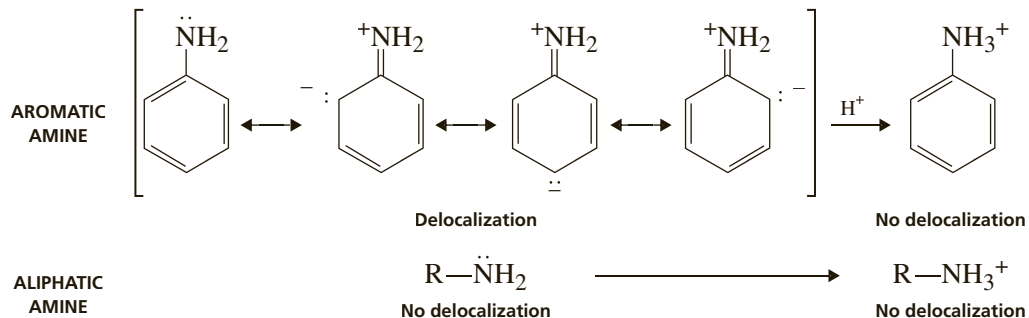





Solubility in 5% NaHCO_3 and 5% NaOH

Compounds that dissolve in sodium bicarbonate, a weak base, are strong acids. Compounds that dissolve in sodium hydroxide, a strong base, may be either strong or weak acids. Thus, one can distinguish weak and strong acids by determining



their solubility in both strong (NaOH) and weak (NaHCO₃) base. The classification of some functional groups as either weak or strong acids is given in the table below.

In this experiment, carboxylic acids (pK_a ~ 5) are generally indicated when a compound is soluble in both bases, and phenols (pK_a ~ 10) are indicated when it is soluble in NaOH only.

Compounds dissolve in base because they form sodium salts that are soluble in the aqueous medium. The salts of some high-molecular-weight compounds are not soluble, however, and precipitate. The salts of the long-chain carboxylic acids, such as myristic acid C₁₄, palmitic acid C₁₆, and stearic acid C₁₈, which form soaps, belong to this category. Some phenols also produce insoluble sodium salts, and often these are colored due to resonance in the anion.

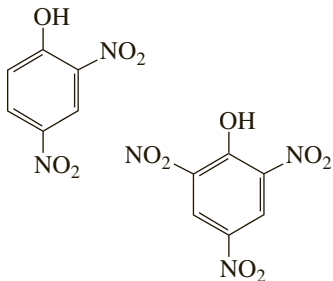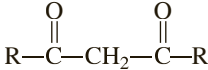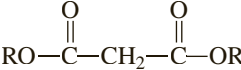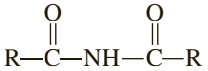Strong Acids (Soluble in Both NaOH and NaHCO ₃)	Weak Acids (Soluble in NaOH but Not in NaHCO ₃)
Sulfonic acids RSO ₃ H Carboxylic acids RCOOH <i>ortho</i> - and <i>para</i> -substituted di- and trinitrophenols 	Phenols ArOH Nitroalkanes RCH ₂ NO ₂ R ₂ CHNO ₂ β-Diketones  β-Diesters  Imides  Sulfonamides ArSO ₂ NH ₂ ArSO ₂ NHR

Both phenols and carboxylic acids produce resonance-stabilized conjugate bases. Thus, bases of appropriate strength may easily remove their acidic protons to form the sodium salts.

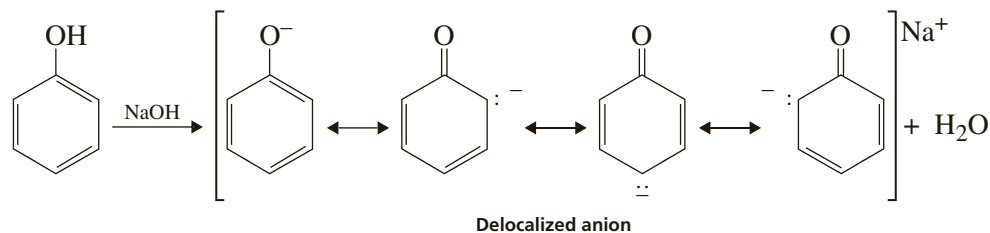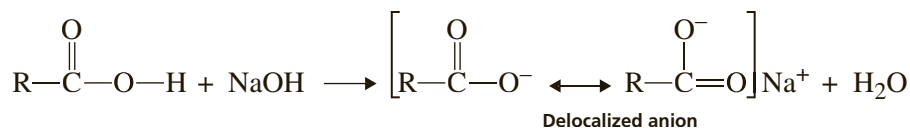


In phenols, substitution of nitro groups in the *ortho* and *para* positions of the ring increases acidity. Nitro groups in these positions provide additional delocalization in the conjugate anion. Phenols that have two or three nitro groups in the *ortho* and *para* positions often dissolve in *both* sodium hydroxide and sodium bicarbonate solutions.

Solubility in Concentrated Sulfuric Acid

Most compounds are soluble in cold, concentrated sulfuric acid, except for inert compounds (alkanes, alkyl halides, and aromatic compounds without a functional group). Of the compounds included in this experiment, alcohols, ketones, aldehydes, and esters belong to this category of being soluble in sulfuric acid. These compounds are described as being “neutral.” However, it must be mentioned that all of the classes of compounds shown in the Solubility chart, except for the inert compounds, will react with cold, concentrated sulfuric acid—not just the “neutral” compounds, alcohols, ketones, aldehydes, and esters. *You must be certain that you have first screened your unknown for solubility in water, sodium bicarbonate, sodium hydroxide, and hydrochloric acid, before trying concentrated sulfuric acid. Inert compounds, which are not included as unknowns, are not soluble in sulfuric acid. Sulfuric acid should always be your last choice of solubility reagents to try.*

Other compounds, not included in this experiment, that also dissolve in sulfuric acid include alkenes, alkynes, ethers, nitroaromatics, and amides. Compounds that are soluble in concentrated sulfuric acid but not in dilute acid are extremely weak bases. Almost any compound containing a nitrogen, an oxygen, or a sulfur atom can be protonated in concentrated sulfuric acid. The ions produced are soluble in the medium.

The classical definition of solubility would be for you to observe the presence of only *one phase* when one of these compounds is mixed with sulfuric acid. However, not all of these compounds will fit this definition. You should expand the definition of “solubility” to include a color change (yellow, brown, or black), the formation of two phases, or possibly even the formation of a precipitate. Any of these observations may be interpreted as “dissolving” in concentrated sulfuric acid. In other words, you should observe some sort of change in the original appearance of the compound.

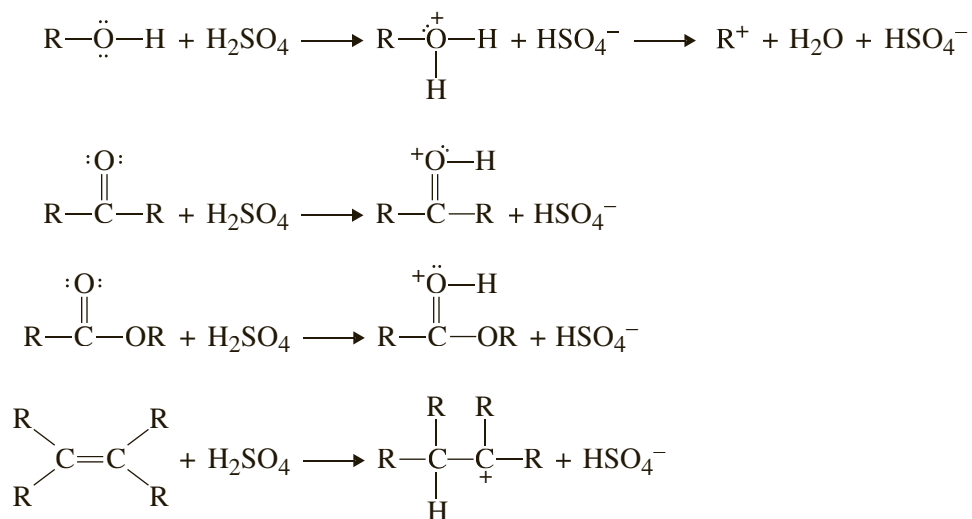


Inert Compounds

Compounds not soluble in concentrated sulfuric acid or any of the other solvents are said to be **inert**. Compounds not soluble in concentrated sulfuric acid include the alkanes, the most simple aromatics, and the alkyl halides. Some examples of inert compounds are hexane, benzene, chlorobenzene, chlorohexane, and toluene.

57B EXPERIMENT 57B

Tests for the Elements (N, S, X)

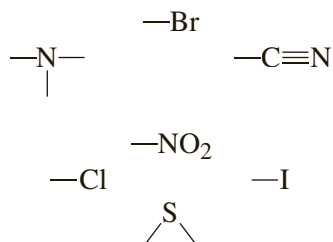

Except for amines (Experiment 57G), which are easily detected by their solubility behavior, all compounds issued in this experiment will contain heteroelements (N, S, Cl, Br, or I) only as *secondary* functional groups. These will be subsidiary to some other important functional group. Thus, no alkyl or aryl halides, nitro compounds, thiols, or thioethers will be issued. However, some of the unknowns may contain a halogen or a nitro group. Less frequently, they may contain a sulfur atom or a cyano group.

Consider as an example *p*-bromobenzaldehyde, an **aldehyde** that contains bromine as a ring substituent. The identification of this compound would hinge on whether the investigator could identify it as an aldehyde. It could probably be identified *without* proving the existence of bromine in the molecule. That information, however, could make the identification easier. In this experiment, methods are given for identifying the presence of a halogen or a nitro group in an unknown compound. Also given is a general method (sodium fusion) for detecting the principal heteroelements that may exist in organic molecules.

Classification tests

<i>Halides</i>	<i>Nitro Groups</i>	<i>N, S, X (Cl, Br, I)</i>
Beilstein test Silver nitrate Sodium iodide/acetone	Ferrous hydroxide	Sodium fusion

SUGGESTED WASTE DISPOSAL

Dispose of all solutions containing silver into a waste container designated for this purpose. Any other aqueous solutions should be disposed of in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container under the hood. This is particularly true of any solution containing benzyl bromide, which is a lachrymator.

TESTS FOR A HALIDE

Beilstein Test

Procedure

Adjust the air and gas mixture so that the flame of a Bunsen burner or microburner is blue. Bend the end of a piece of copper wire so that a small closed loop is created.

Heat the loop end of the wire in the flame until it glows brightly. After the wire has cooled, dip the wire directly into a sample of the unknown. If the unknown is a solid and won't adhere to the copper wire, place a small amount of the substance on a watch glass, wet the copper wire in distilled water, and place the wire into the sample on the watch glass. The solid should adhere to the wire. Again heat the wire in the Bunsen burner flame. The compound will first burn. After the burning, a green flame will be produced if a halogen is present. You should hold the wire in the flame either just above the tip of the flame or at its outside edge near the bottom of the flame. You will need to experiment to find the best position to hold the copper wire to obtain the best result.

Test Compounds

Try this test on bromobenzene and benzoic acid.

Discussion

Halogens can be detected easily and reliably by the Beilstein test. It is the simplest method for determining the presence of a halogen, but it does not differentiate among chlorine, bromine, and iodine, any one of which will give a positive test. However, when the identity of the unknown has been narrowed to two choices, of which one has a halogen and one does not, the Beilstein test will often be enough to distinguish between the two.

A positive Beilstein test results from the production of a volatile copper halide when an organic halide is heated with copper oxide. The copper halide imparts a blue-green color to the flame.

This test can be very sensitive to small amounts of halide impurities in some compounds. Therefore, use caution in interpreting the results of the test if you obtain only a weak color.

Silver Nitrate Test

Procedure

Add 1 drop of a liquid or 5 drops of a concentrated ethanolic solution of the unknown solid to 2 mL of a 2% ethanolic silver nitrate solution. If no reaction is observed after 5 minutes at room temperature, heat the solution in a hot water bath at about 100°C, and note whether a precipitate forms. If a precipitate forms, add 2 drops of 5% nitric acid, and note whether the precipitate dissolves. Carboxylic acids give a false test by precipitating in silver nitrate, but they dissolve when nitric acid is added. Silver halides, in contrast, do not dissolve in nitric acid.

Test Compounds

Apply this test to benzyl bromide (α -bromotoluene) and bromobenzene. Discard all waste reagents in a suitable waste container in the hood because benzyl bromide is a lachrymator.


Discussion

This test depends on the formation of a white or off-white precipitate of silver halide when silver nitrate is allowed to react with a sufficiently reactive halide.



The test does not distinguish among chlorides, bromides, and iodides but does distinguish labile (reactive) halides from halides that are unreactive. Halides substituted on an aromatic ring will not usually give a positive silver nitrate test; however, alkyl halides of many types will give a positive test.

The most reactive compounds are those able to form stable carbocations in solution and those equipped with good leaving groups ($X = \text{I}, \text{Br}, \text{Cl}$). Benzyl, allyl, and tertiary halides react immediately with silver nitrate. Secondary and primary halides do not react at room temperature but react readily when heated. Aryl and vinyl halides do not react at all, even at elevated temperatures. This pattern of reactivity fits the stability order for various carbocations quite well. Compounds that produce stable carbocations react at higher rates than those that do not.

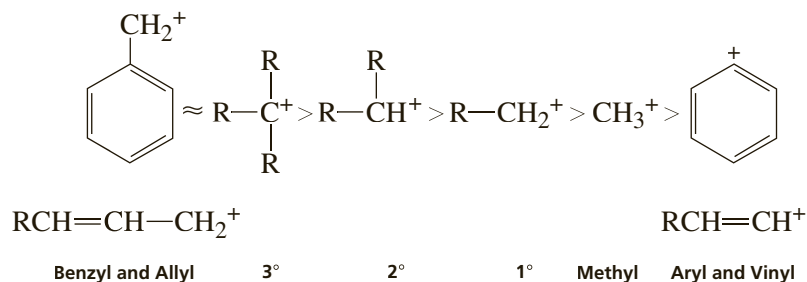

The fast reaction of benzylic and allylic halides is a result of the resonance stabilization that is available to the intermediate carbocations formed. Tertiary halides are more reactive than secondary halides, which are in turn more reactive than primary or methyl halides because alkyl substituents are able to stabilize the intermediate carbocations by an electron-releasing effect. The methyl carbocations have no alkyl groups and are the least stable of all carbocations mentioned thus far. Vinyl and aryl carbocations are extremely unstable because the charge is localized on an sp^2 -hybridized carbon (double-bond carbon) rather than one that is sp^3 -hybridized.

Sodium Iodide In Acetone

Procedure

Add 4 drops of the unknown to 2 mL of a 15% NaI-in-acetone solution in a small test tube. Do not use your thumb or a stopper, but mix the solution in the following manner. Hold the top of the test tube between the thumb and index finger of one hand and “flick” the bottom of the test tube using the index finger of your other hand. Repeat this several times to assure adequate mixing. If a precipitate appears, note the amount of time required for it to appear.

Repeat the mixing from time to time over a period of about 5 minutes. Any color changes should be ignored. Appearance of a precipitate is the significant result.

If no precipitate appears after the 5-minute period, heat the test tube in a 50°C water bath for a few minutes. Do not exceed 50°C, or the acetone solvent will boil away.

Discard all waste reagents in a suitable waste container in the hood because many organic halides are lachrymators.

Test Compounds

Apply this test to benzyl bromide (α -bromotoluene) and bromobenzene.

Discussion

Generally reactive halides give a precipitate within 3 minutes at room temperature, moderately reactive halides give a precipitate when heated, and unreactive halides do not give a precipitate, even after being heated.

A reagent composed of sodium or potassium iodide dissolved in acetone is useful in classifying alkyl halides according to their reactivity in an S_N2 reaction.

Iodide is an excellent nucleophile, and acetone is a nonpolar solvent. The tendency to form a precipitate increases the completeness of the reaction. Sodium iodide and potassium iodide are soluble in acetone, but sodium and potassium chlorides and bromides are not. Consequently, as bromide ion or chloride ion is produced, the sodium salt is precipitated from the solution. According to LeChâtelier's Principle, the precipitation of a product from the reaction solution drives the equilibrium toward the right; such is the case in the reaction described here:



Recall that allylic, benzylic, and primary halides react readily in $\text{S}_{\text{N}}2$ reactions, secondary halides are slower, while tertiary and sterically hindered halides react very slowly or not at all. Vinyl and aromatic halides generally do not react.

DETECTION OF NITRO GROUPS

Although nitro compounds will not be issued as distinct unknowns, many of the unknowns may have a nitro group as a secondary functional group. The presence of a nitro group, and hence nitrogen, in an unknown compound is determined most easily by infrared spectroscopy. Unfortunately, functional groups other than the nitro group may also give a positive result. You should interpret the results of this test with caution. It is very important to perform this procedure on the test compound. It may also help to do this on a known compound without a nitro group.

Ferrous Hydroxide Test

Procedure


Place 1.5 mL of freshly prepared 5% aqueous ferrous ammonium sulfate in a small test tube, and add about 10 mg of a solid or 5 drops of a liquid compound. Mix the solution well, and then add first 1 drop of 2 M sulfuric acid and then 1 mL of 2 M potassium hydroxide in methanol. Stopper the test tube and shake it vigorously. A positive test is indicated by the formation of a red-brown precipitate, usually within 1 minute.

Test Compound

Apply this test to 2-nitrotoluene.

Discussion

Most nitro compounds oxidize ferrous hydroxide to ferric hydroxide, which is a red-brown solid. A precipitate indicates a positive test.

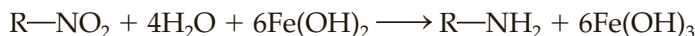


Infrared Spectroscopy

The nitro group gives two strong bands near 1560 cm^{-1} and 1350 cm^{-1} . See Technique 25 for details.

DETECTION OF A CYANO GROUP

Although nitriles will not be given as unknowns in this experiment, the cyano group may be a subsidiary functional group whose presence or absence is important to the final identification of an unknown compound. The cyano group can

be hydrolyzed in a strong base by heating vigorously to give carboxylic acid and ammonia gas:

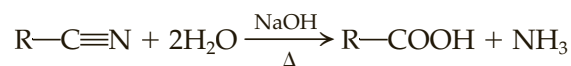

The ammonia gas can be detected by its odor or by using moist pH paper. However, this method is somewhat difficult, and the presence of a nitrile group is confirmed most easily by infrared spectroscopy. No other functional groups (except some $\text{C}\equiv\text{C}$) absorb in the same region of the spectrum as $\text{C}\equiv\text{N}$.

Infrared Spectroscopy

$\text{C}\equiv\text{N}$ stretch is a sharp band of medium intensity near 2250 cm^{-1} . See Technique 25 for details.

SODIUM FUSION TEST (OPTIONAL)

When an organic compound containing nitrogen, sulfur, or halide atoms is fused with sodium metal, there is a reductive decomposition of the compound, which converts these atoms to the sodium salts of the inorganic ions CN^- , S^{2-} , and X^- .



When the fusion mixture is dissolved in distilled water, the cyanide, sulfide, and halide ions can be detected by standard qualitative inorganic tests.

CAUTION



Always remember to manipulate the sodium metal with a knife or a forceps. Do not touch it with your fingers. Keep sodium away from water. Destroy all waste sodium with 1-butanol or ethanol. Wear safety glasses.

General Method for Preparing Stock Solution

Procedure

Using a forceps and a knife, take some sodium from the storage container, cut a small piece about the size of a small pea (3 mm on a side), and dry it on a paper towel. Place this small piece of sodium in a clean, dry, small test tube (10 mm × 75 mm). Clamp the test tube to a ring stand, and heat the bottom of the tube with a microburner until the sodium melts and its metallic vapor can be seen to rise about a third of the way up the tube. The bottom of the tube will probably have a dull red glow. Remove the burner and *immediately* drop the sample directly into the tube. Use about 10 mg of a solid placed on the end of a spatula or 2–3 drops of a liquid. Be sure to drop the sample directly down the center of the tube so that it touches the hot sodium metal and does not adhere to the side of the test tube. If the fusion is successful, there will usually be a flash or a small explosion. If the reaction is not successful, heat the tube to red heat for a few seconds to ensure complete reaction.

Allow the test tube to cool to room temperature, and then carefully add 10 drops of methanol, a drop at a time, to the fusion mixture. Using a spatula or a long glass rod, reach into the test tube and stir the mixture to ensure complete reaction of any excess sodium metal. The fusion will have destroyed the test tube for other uses. Thus, the easiest way to recover the fusion mixture is to crush the test tube into a small beaker containing 5–10 mL of *distilled* water. The tube is easily crushed if it is

placed in the angle of a clamp holder. Tighten the clamp until the tube is securely held near its bottom and then—standing back from the beaker and holding the clamp at its opposite end—continue tightening the clamp until the test tube breaks and the pieces fall into the beaker. Stir the solution well, heat until it boils, and then filter it by gravity through a fluted filter (see Technique 8, Figure 8.3). Portions of this solution will be used in the tests to detect nitrogen, sulfur, and the halogens.

This stock solution or the one prepared by the alternate method below will be used in the nitrogen, sulfur, and halides tests that follow.

Alternative Method for Preparing Stock Solution

Procedure

With some volatile liquids, the previous method will not work. The compounds volatilize before they reach the sodium vapors. For such compounds, place 4 or 5 drops of the pure liquid in a clean, dry test tube, clamp it, and cautiously add the small piece of sodium metal. If there is any reaction, wait until it subsides. Then heat the test tube to red heat, and continue according to the instructions in the second paragraph of the preceding procedure. This stock solution or the one prepared by the general method above will be used in the nitrogen, sulfur, and halides tests that follow.

Nitrogen Test

Procedure

Using pH paper and a 10% sodium hydroxide solution, adjust the pH of about 1 mL of the stock solution to pH 13. Add 2 drops of saturated ferrous ammonium sulfate solution and 2 drops of 30% potassium fluoride solution. Boil the solution for about 30 seconds. Then acidify the hot solution by adding 30% sulfuric acid dropwise until the iron hydroxides dissolve. Avoid using excess acid. If nitrogen is present, a dark Prussian blue (not green) precipitate $\text{NaFe}_2(\text{CN})_6$ will form, or the solution will assume a dark blue color.

Reagents

Dissolve 5 g of ferrous ammonium sulfate in 100 mL of water. Dissolve 30 g of potassium fluoride in 100 mL of water.

Sulfur Test

Procedure

Acidify about 1 mL of the stock solution prepared above with acetic acid, and add a few drops of a 1% lead acetate solution. The presence of sulfur is indicated by a black precipitate of lead sulfide (PbS).

CAUTION



Many compounds of lead(II) are suspected carcinogens (see Technique 1, Section 1.4) and should be handled with care. Avoid contact.


Halide Tests

Procedure

Cyanide and sulfide ions interfere with the test for halides. If such ions are present, they must be removed. To accomplish this, acidify 2 mL of the stock solution prepared above with dilute nitric acid and boil it for about 2 minutes. This will drive off any HCN or H_2S that is formed. When the solution cools, add a few drops of a 5% silver nitrate solution. A *voluminous* precipitate indicates a halide. A faint turbidity *does not* mean a positive test. Silver chloride is white. Silver bromide is off-white. Silver iodide is yellow. Silver chloride will readily dissolve in concentrated ammonium hydroxide, whereas silver bromide is only slightly soluble.

Differentiation of Chloride, Bromide, and Iodide**Procedure**

Acidify 2 mL of the stock solution prepared above with 10% sulfuric acid, and boil it for about 2 minutes. Cool the solution and add about 0.5 mL of methylene chloride. Add a few drops of chlorine water or 2–4 mg of calcium hypochlorite.¹ Check to be sure that the solution is still acidic. Then stopper the tube, shake it vigorously, and set it aside to allow the layers to separate. An orange to brown color in the methylene chloride layer indicates bromine. Violet indicates iodine. No color or a light yellow indicates chlorine.

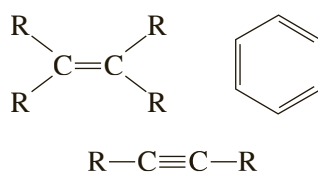57C EXPERIMENT 57C*Tests for Unsaturation*

The unknowns to be issued for this experiment have neither a double bond nor a triple bond as their *only* functional group. Hence, simple alkenes and alkynes can be ruled out as possible compounds. Some of the unknowns may have a double or a triple bond, however, *in addition to* another more important functional group. The tests described allow you to determine the presence of a double bond or a triple bond (unsaturation) in such compounds.

Classification tests

<i>Unsaturation</i>	<i>Aromaticity</i>
Bromine–methylene chloride	Ignition test
Potassium permanganate	

SUGGESTED WASTE DISPOSAL

Test reagents that contain bromine should be discarded into a special waste container designated for this purpose. Methylene chloride must be placed in the organic waste container designated for the disposal of halogenated organic wastes. Dispose of all other aqueous solutions in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container.

¹ Clorox, the commercial bleach, is a permissible substitute for chlorine water, as is any other brand of bleach, provided that it is based on sodium hypochlorite.

TEST FOR SIMPLE MULTIPLE BONDS

Bromine in Methylene Chloride

Procedure

Dissolve 50 mg of the unknown solid or 4 drops of the unknown liquid in 1 mL of methylene chloride (dichloromethane) or in 1,2-dimethoxyethane. Add a 2% (by volume) solution of bromine in methylene chloride, dropwise, with shaking. If you find that the red color remains after adding 1 or 2 drops of the bromine solution, the test is negative. If the red color disappears, continue adding the bromine in methylene chloride until the red bromine color remains. The test is positive if more than 5 drops of the bromine solution were added, with discharge of the red color of bromine. If the red color disappears, try adding more drops of the bromine solution to see how many drops are necessary before the red color persists. Usually, many drops of the bromine solution will be decolorized when an isolated double bond is present. Hydrogen bromide should not be evolved. If hydrogen bromide gas is evolved, you will note a "fog" when you blow across the mouth of the test tube. The HBr can also be detected by a moistened piece of litmus or pH paper. If hydrogen bromide is evolved, the reaction is a **substitution reaction** (see following discussion) and not an **addition reaction**, and a double or triple bond is probably not present.

Reagent


The classic method for running this test is to use bromine dissolved in carbon tetrachloride. Because of the toxic nature of this solvent, methylene chloride has been substituted for carbon tetrachloride. The instructor must prepare this reagent because of the danger associated with the very toxic bromine vapor. Be sure to work in an efficient fume hood. Dissolve 2 mL of bromine in 100 mL of methylene chloride (dichloromethane). The solvent will undergo a light-induced, free-radical substitution producing hydrogen bromide over a period of time. After about 1 week, the color of the 2% solution of bromine in methylene chloride fades noticeably, and the odor of the HBr can be detected in the reagent. Although the decolorization tests still work satisfactorily, the presence of HBr makes it difficult to distinguish between addition and substitution reactions. A freshly prepared solution of bromine in methylene chloride must be used to make this distinction. Deterioration of the reagent can be forestalled by storing it in a brown glass bottle.

Test Compounds

Try this test with cyclohexene, cyclohexane, toluene, and acetone.


Discussion

A successful test depends on the addition of bromine, a red liquid, to a double or a triple bond to give a colorless dibromide:

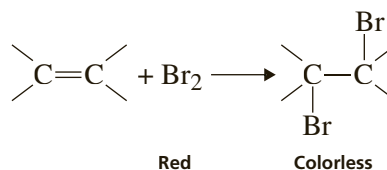


Not all double bonds react with the bromine solution. Only those that are electron-rich are sufficiently reactive nucleophiles to initiate the reaction. A double bond

that is substituted by electron-withdrawing groups often fails to react or reacts slowly. Fumaric acid is an example of a compound that fails to give the reaction.

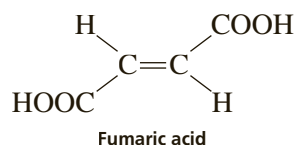

Aromatic compounds either do not react with the bromine reagent, or they react by **substitution**. Only the aromatic rings that have activating groups as substituents ($-\text{OH}$, $-\text{OR}$, or $-\text{NR}_2$) give the substitution reaction.

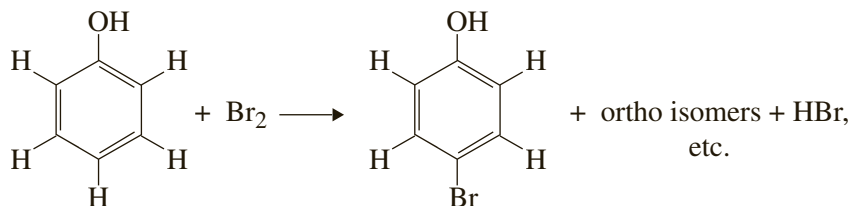

Some ketones and aldehydes react with bromine to give a **substitution product**, but this reaction is slow except for ketones that have a high enol content. When substitution occurs, not only is the bromine color discharged, but hydrogen bromide gas is also evolved.

Potassium Permanganate (Baeyer Test)

Procedure


Dissolve 25 mg of the unknown solid or 2 drops of the unknown liquid in 2 mL of 95% ethanol (1,2-dimethoxyethane may also be used). Slowly add a 1% aqueous solution (weight/volume) of potassium permanganate, drop by drop while shaking, to the unknown. In a positive test, the purple color of the reagent is discharged, and a brown precipitate of manganese dioxide forms, usually within 1 minute. If alcohol was the solvent, the solution should not be allowed to stand for more than 5 minutes, because oxidation of the alcohol will begin slowly. Because permanganate solutions undergo some decomposition to manganese dioxide on standing, any small amount of precipitate should be interpreted with caution.

Test Compounds

Try this test on cyclohexene and toluene.

Discussion

This test is positive for double and triple bonds but not for aromatic rings. It depends on the conversion of the purple ion MnO_4^- to a brown precipitate of MnO_2 following the oxidation of an unsaturated compound.

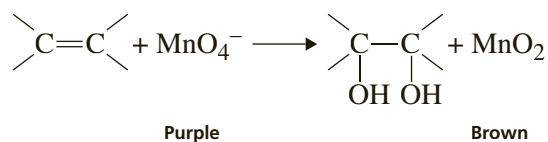

Other easily oxidized compounds also give a positive test with potassium permanganate solution. These substances include aldehydes, some alcohols, phenols, and aromatic amines. If you suspect that any of these functional groups is present, you should interpret the test with caution.

Spectroscopy*Infrared***Double Bonds (C=C)**

C=C stretch usually occurs near 1680–1620 cm^{-1} .

Symmetrical alkenes may have no absorption.

C—H stretch of vinyl hydrogens occurs $> 3000 \text{ cm}^{-1}$, but usually not higher than 3150 cm^{-1} .

C—H out-of-plane bending occurs near 1000–700 cm^{-1} .

Triple Bonds (C≡C)

C≡C stretch usually occurs near 2250–2100 cm^{-1} . The peak is usually sharp. Symmetrical alkynes show no absorption.

C—H stretch of terminal acetylenes occurs near 3310–3200 cm^{-1} .

See Technique 25 for details.

Nuclear Magnetic Resonance

Vinyl hydrogens have resonance near 5–7 ppm and have coupling values as follows: $J_{trans} = 11\text{--}18 \text{ Hz}$, $J_{cis} = 6\text{--}15 \text{ Hz}$, $J_{geminal} = 0\text{--}5 \text{ Hz}$. Allylic hydrogens have resonance near 2 ppm. Acetylenic hydrogens have resonance near 2.8–3.0 ppm. See Technique 26 for details on proton NMR. Carbon NMR is described in Technique 27.

TESTS FOR AROMATICITY

None of the unknowns to be issued for this experiment will be simple aromatic hydrocarbons. All aromatic compounds will have a principal functional group as a part of their structure. Nevertheless, in many cases it will be useful to be able to recognize the presence of an aromatic ring. Although infrared and nuclear magnetic spectroscopy provide the most reliable methods of determining aromatic compounds, often they can be detected by a simple ignition test.

Ignition Test*Procedure*

Working in a hood, place a small amount of the compound on a spatula and place it in the flame of a Bunsen burner. Observe whether a sooty flame results. Compounds giving the sooty yellow flame have a high degree of unsaturation and may be aromatic. This test should be interpreted with care because some nonaromatic compounds may produce soot. If in doubt, use spectroscopy to more reliably determine the presence or absence of an aromatic ring.

Test Compounds

Try this test with ethyl benzoate and benzoin.

Discussion

The presence of an aromatic ring will usually lead to the production of a sooty yellow flame in this test. In addition, halogenated alkanes and high-molecular-weight aliphatic compounds may produce a sooty yellow flame. Aromatic compounds with high oxygen content may burn cleaner and produce less soot even though they contain an aromatic ring.

This is actually a test to determine the ratio of carbon to hydrogen, and oxygen in an unknown substance. If the carbon-to-hydrogen ratio is high and if little or no oxygen is present, you will observe a sooty flame. For instance, acetylene, C_2H_2 (a gas), will burn with a sooty flame unless mixed with oxygen. When the carbon-to-hydrogen ratio is nearly equal to one, you will be very likely to see a sooty flame.

Spectroscopy**Infrared**

C=C aromatic-ring double bonds appear in the 1600–1450 cm^{-1} region. There are often four sharp absorptions that occur in pairs near 1600 cm^{-1} and 1450 cm^{-1} , which are characteristic of an aromatic ring.

Special ring absorptions: There are often weak ring absorptions around 2000–1600 cm^{-1} . These are frequently obscured, but when they can be observed, the relative shapes and numbers of these peaks can often be used to ascertain the type of ring substitution.


=C—H stretch, aromatic ring: The aromatic C—H stretch always occurs at a higher frequency than 3000 cm^{-1} .

=C—H out-of-plane bending peaks appear in the region 900–690 cm^{-1} . The number and position of these peaks can be used to determine the substitution pattern of the ring.

See Technique 25 for details.

Nuclear Magnetic Resonance

Hydrogens attached to an aromatic ring usually have resonance near 7 ppm. Mono-substituted rings not substituted by anisotropic or electronegative groups often give a single resonance for all of the ring hydrogens. Monosubstituted rings with anisotropic or electronegative groups usually have the aromatic resonances split into two groups integrating either 3:2 or 2:3. A nonsymmetric, *para*-disubstituted ring has a characteristic four-peak splitting pattern (see Technique 26). Carbon NMR is described in Technique 27.

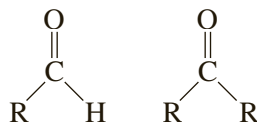57D EXPERIMENT 57D**Aldehydes and Ketones**

Compounds containing the carbonyl functional group $\text{C}=\text{O}$, where it has only hydrogen atoms or alkyl groups as substituents, are called aldehydes, RCHO , or ketones, RCOR' . The chemistry of these compounds is primarily due to the chemistry of the carbonyl functional groups. These compounds are identified by the distinctive reactions of the carbonyl function.

Solubility Characteristics

HCl	NaHCO_3	NaOH	H_2SO_4	Ether
(-)	(-)	(-)	(+)	(+)
Water: < C_5 and some C_6 (+)				
> C_5 (-)				

Classification Tests

Aldehydes and ketones	
2,4-Dinitrophenylhydrazine	
Aldehydes only	Methyl ketones
Tollens reagent	Iodoform test
Chromic acid	
Compounds with high enol content	
Ferric chloride test	

SUGGESTED WASTE DISPOSAL

Solutions containing 2,4-dinitrophenylhydrazine or derivatives formed from it should be placed in a waste container designated for these compounds. Any solution containing chromium must be disposed of in a waste container specifically identified for the disposal of chromium wastes. Dispose of all solutions containing silver by acidifying them with 5% hydrochloric acid and then placing them in a waste container designated for this purpose. Dispose of all other aqueous solutions in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container.

CLASSIFICATION TESTS

Most aldehydes and ketones give a solid, yellow to red precipitate when mixed with 2,4-dinitrophenylhydrazine. However, only aldehydes will reduce chromium(VI) or silver(I). By this difference in behavior, you can differentiate between aldehydes and ketones.

2,4-Dinitrophenylhydrazine

Procedure

Place 1 drop of the liquid unknown in a small test tube and add 1 mL of the 2,4-dinitrophenylhydrazine reagent. If the unknown is a solid, dissolve about 10 mg (estimate) in a minimum amount of 95% ethanol or di(ethylene glycol) diethyl ether before adding the reagent. Shake the mixture vigorously. Most aldehydes and ketones will give a yellow to red precipitate immediately. However, some compounds will require up to 15 minutes, or even *gentle* heating, to give a precipitate. A precipitate indicates a positive test.

Test Compounds

Try this test on cyclohexanone, benzaldehyde, and benzophenone.

CAUTION



Many derivatives of phenylhydrazine are suspected carcinogens (see Technique 1, Section 1.4) and should be handled with care. Avoid contact.


Reagent

Dissolve 3.0 g of 2,4-dinitrophenylhydrazine in 15 mL of concentrated sulfuric acid. In a beaker, slowly add, with mixing, 23 mL of water until the solid dissolves. Add 75 mL of 95% ethanol to the warm solution, while stirring. After thorough mixing, filter the solution if any solid remains. This reagent needs to be prepared fresh each time.

Discussion

Most aldehydes and ketones give a precipitate, but esters generally do not give this result. Thus, an ester usually can be eliminated by this test. The color of the 2,4-dinitrophenylhydrazone (precipitate) formed is often a guide to the amount of conjugation in the original aldehyde or ketone. Unconjugated ketones, such as cyclohexanone, give yellow precipitates, whereas conjugated ketones, such as benzophenone, give orange to red precipitates. Compounds that are highly conjugated give red precipitates. However, the 2,4-dinitrophenylhydrazine reagent is itself orange-red, and the color of any precipitate must be judged cautiously.

Occasionally, compounds that are either strongly basic or strongly acidic precipitate the unreacted reagent.

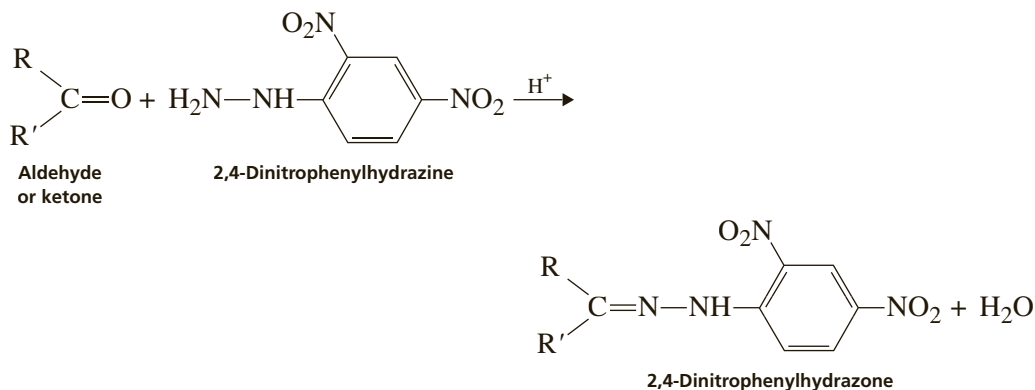

Some allylic and benzylic alcohols give this test result because the reagent can oxidize them to aldehydes and ketones, which subsequently react. Some alcohols may be contaminated with carbonyl impurities, either as a result of their method of synthesis (reduction) or as a result of their becoming air-oxidized. A precipitate formed from small amounts of impurity in the solution will be formed in small amounts. With some caution, a test that gives only a slight amount of precipitate can usually be ignored. The infrared spectrum of the compound should establish its identity and identify any impurities present.

Tollens Test

Procedure

The reagent must be prepared immediately before use. To prepare the reagent, mix 1 mL of Tollens solution A with 1 mL of Tollens solution B. A precipitate of silver oxide will form. Add enough dilute (10%) ammonia solution (dropwise) to the mixture to dissolve the silver oxide *just barely*. The reagent so prepared can be used immediately for the following test.

Dissolve 1 drop of a liquid aldehyde or 10 mg (approximate) of a solid aldehyde in the minimum amount of di(ethylene glycol) diethyl ether. Add this solution, a little at a time, to the 2–3 mL of reagent contained in a small test tube. Shake the solution well. If a mirror of silver is deposited on the inner walls of the test tube, the test is positive. In some cases, it may be necessary to warm the test tube in a warm-water bath.

Test Compounds

Try the test on benzaldehyde, butanal (butyraldehyde), and cyclohexanone.

CAUTION



The reagent should be prepared immediately before use and all residues disposed of immediately after use. Dispose of any residues by acidifying them with 5% hydrochloric acid and then placing them in a waste container designated for this purpose. On standing, the reagent tends to form silver fulminate, a *very explosive* substance. Solutions containing the mixed Tollens reagent should never be stored.

Reagents

Solution A: Dissolve 3.0 g of silver nitrate in 30 mL of water. *Solution B:* Prepare a 10% sodium hydroxide solution.

Discussion

Most aldehydes reduce ammoniacal silver nitrate solution to give a precipitate of silver metal. The aldehyde is oxidized to a carboxylic acid:



Ordinary ketones do not give a positive result in this test. The test should be used only if it has already been shown that the unknown compound is either an aldehyde or a ketone.

Chromic Acid Test: Alternative Test

CAUTION



Many chromium (VI) compounds are suspected carcinogens. If you would like to run this test, talk to your instructor first. Most often, the Tollens test will easily distinguish between aldehydes and ketones, and you should do that test first. If you run the chromic acid test, be sure to wear gloves to avoid contact with this reagent.

Procedure

Dissolve 1 drop of a liquid or 10 mg (approximate) of a solid aldehyde in 1 mL of *reagent-grade* acetone. Add several drops of the chromic acid reagent, a drop at a time, while shaking the mixture. A positive test is indicated by a green precipitate and a loss of the orange color in the reagent. With aliphatic aldehydes, RCHO, the solution turns cloudy within 5 seconds, and a precipitate appears within 30 seconds. With aromatic aldehydes, ArCHO, it generally takes 30–120 seconds for a precipitate to form, but with some it may take even longer. In some cases, however, you may find that some of the original orange color may remain, together with a green or brown precipitate. This should be interpreted as a positive test. In a negative test, a nongreen precipitate may form in an orange solution.

In performing this test, make sure that the acetone used for the solvent does not give a positive test with the reagent. Add several drops of the chromic acid reagent to a few drops of the reagent acetone contained in a small test tube. Allow this mixture to stand for 3–5 minutes. If no reaction has occurred by this time, the acetone is pure enough to use as a solvent for the test. If a positive test resulted, try another bottle of acetone.

Test Compounds


Try the test on benzaldehyde, butanal (butyraldehyde), and cyclohexanone.

Reagent

Dissolve 20 g of chromium trioxide (CrO_3) in 60 mL of cold water in a beaker. With stirring, slowly and carefully add 20 mL of concentrated sulfuric acid to the solution. This reagent should be prepared fresh each time.

Discussion

This test is based on the fact that aldehydes are easily oxidized to the corresponding carboxylic acid by chromic acid. The green precipitate is due to chromous sulfate.

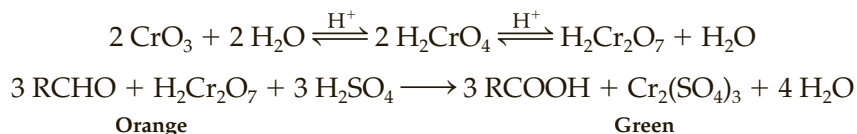

Primary and secondary alcohols are also oxidized by this reagent (see Experiment 57H). Therefore, this test is not useful in identifying aldehydes *unless* a positive identification of the carbonyl group has already been made. Aldehydes give a 2,4-dinitrophenylhydrazine test result, whereas alcohols do not.

There are numerous other tests used to detect the aldehyde functional group. Most are based on an easily detectable oxidation of the aldehyde to a carboxylic acid. The most common tests are the Tollens, Fehling's, and Benedict's tests. Only the Tollens test is described in this book. The Tollens test is often more reliable than the chromic acid test for aldehydes.

Iodoform Test

Procedure

Prepare a 60°C–70°C water bath in a beaker. Using a Pasteur pipette, add 6 drops of a liquid unknown to a 15-mm × 100-mm or 15-mm × 125-mm test tube. Alternatively, 0.06 g of the unknown solid may be used. Dissolve the unknown liquid or solid compound in 2 mL of 1,2-dimethoxyethane. Add 2 mL of 10% aqueous sodium hydroxide solution, and place the test tube in the hot-water bath. Next add 4 mL of iodine–potassium iodide solution in 1-mL portions to the test tube. *Cork* the test tube and shake it after adding each portion of iodine reagent. Heat the mixture in the hot-water bath for about 5 minutes, shaking the test tube occasionally. It is likely that some or all of the dark color of the iodine reagent will be discharged.

If the dark color of the iodine reagent is still apparent following heating, add 10% sodium hydroxide solution until the dark color of the iodine reagent has been discharged. Shake the mixture in the test tube (corked) during the addition of sodium hydroxide. Care need not be taken to avoid adding excess sodium hydroxide.

After the dark iodine color of the solution has been discharged, fill the test tube with water to within 2 cm of the top. *Cork* the test tube and shake it vigorously. Allow the tube to stand for at least 15 minutes at room temperature. The appearance of a pale yellow precipitate of iodoform, CHI_3 , constitutes a positive test, indicating that the unknown is a methyl ketone or a compound that is easily oxidized to a methyl ketone, such as a 2-alkanol. Other ketones will also decolorize the iodine solution, but they will not give a precipitate of iodoform *unless* there is an impurity of a methyl ketone present in the unknown.

The yellow precipitate usually settles out slowly onto the bottom of the test tube. Sometimes, the yellow color of iodoform is masked by a dark substance. If this is the case, *cork* the test tube and shake it vigorously. If the dark color persists, add more sodium hydroxide solution, and shake the test tube again. Then allow the tube to stand for at least 15 minutes. If there is some doubt as to whether the solid is iodoform, collect the precipitate on a Hirsch funnel and dry it. Iodoform melts at 119°C–121°C.

You may find on some occasions that methyl ketone gives only a yellow coloration to the solution rather than a distinct yellow precipitate. You should be cautious about drawing any conclusions from this result. Therefore, you should depend on proton NMR to confirm the presence of a methyl group attached directly to a carbonyl group (singlet at about 2 ppm).

Test Compounds


Try the test on 2-heptanone, 4-heptanone (dipropyl ketone), and 2-pentanol.

Reagents

The iodine reagent is prepared by dissolving 20 g of potassium iodide and 10 g of iodine in 100 mL of water. The aqueous sodium hydroxide solution is prepared by dissolving 10 g of sodium hydroxide in 100 mL of water.

Discussion

The basis of this test is the ability of certain compounds to form a precipitate of iodoform when treated with a basic solution of iodine. Methyl ketones are the most common types of compounds that give a positive result in this test. However, acetaldehyde, CH_3CHO , and alcohols with the hydroxyl group at the 2-position of the chain also give a precipitate of iodoform. 2-Alkanols of the type described are easily oxidized to methyl ketones under the conditions of the reaction. The other product of the reaction, besides iodoform, is the sodium or potassium salt of a carboxylic acid.

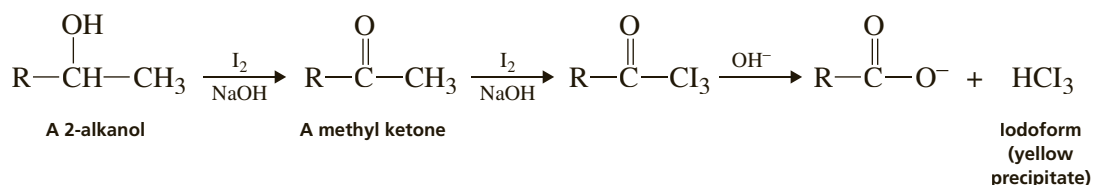

Ferric Chloride Test *Procedure*

Some aldehydes and ketones, those that have a high **enol content**, give a positive ferric chloride test, as described for phenols in Experiment 57F.

Spectroscopy *Infrared*

The carbonyl group is usually one of the strongest-absorbing groups in the infrared spectrum, with a very broad range: $1800\text{--}1650\text{ cm}^{-1}$. The aldehyde functional group has *very characteristic* C—H stretch absorptions: two sharp peaks that lie *far outside* the usual region for —C—H , =C—H or ≡C—H .

Aldehydes

C=O stretch at approximately 1725 cm^{-1} is normal.
 $1725\text{--}1685\text{ cm}^{-1}$.*

C—H stretch (aldehyde—CHO) has two weak bands at about 2750 cm^{-1} and 2850 cm^{-1} .

See Technique 25 for details.

Ketones

C=O stretch at approximately 1715 cm^{-1} is normal.
 $1780\text{--}1665\text{ cm}^{-1}$.*

Nuclear Magnetic Resonance





Hydrogens alpha to a carbonyl group have resonance in the region between 2 ppm and 3 ppm. The hydrogen of an aldehyde group has a characteristic resonance between 9 ppm and 10 ppm. In aldehydes, there is coupling between the aldehyde hydrogen and any alpha hydrogens ($J = 1\text{--}3\text{ Hz}$).

See Technique 26 for details on proton NMR. Carbon NMR is described in Technique 27.

* **Conjugation** moves the absorption to lower frequencies. **Ring strain** (cyclic ketones) moves the absorption to higher frequencies.

Derivatives

The most common derivatives of aldehydes and ketones are 2,4-dinitrophenylhydrazones, oximes, and semicarbazones. Procedures for preparing these derivatives are given in Appendix 2.

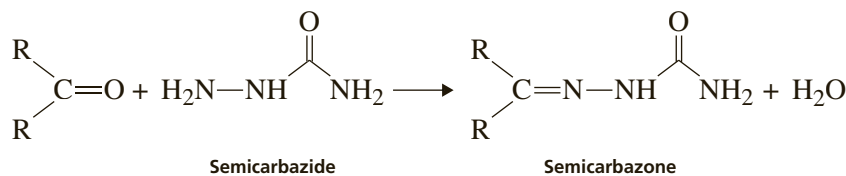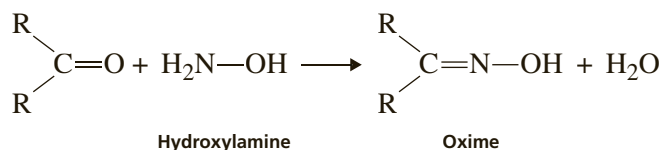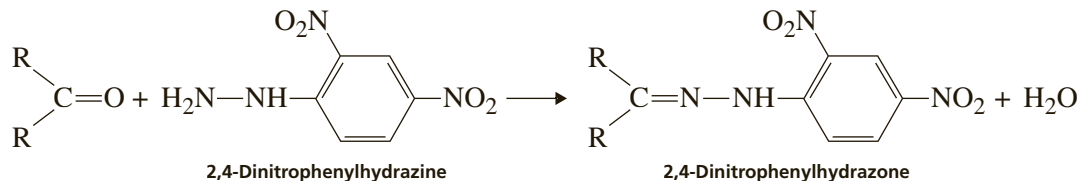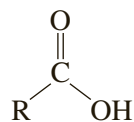**57E****EXPERIMENT 57E****Carboxylic Acids**

Carboxylic acids are detectable mainly by their solubility characteristics. They are soluble in *both* dilute sodium hydroxide and sodium bicarbonate solutions.

Solubility Characteristics					Classification Tests
HCl	NaHCO ₃	NaOH	H ₂ SO ₄	Ether	pH of an aqueous solution
(-)	(+)	(+)	(+)	(+)	Sodium bicarbonate
Water: < C ₆ (+)					Silver nitrate
> C ₆ (-)					Neutralization equivalent

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous solutions in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container.

CLASSIFICATION TESTS

pH of an Aqueous Solution

Procedure


If the compound is soluble in water, simply prepare an aqueous solution and check the pH with pH paper. If the compound is an acid, the solution will have a low pH.

Compounds that are insoluble in water can be dissolved in ethanol (or methanol) and water. First, dissolve the compound in the alcohol, and then add water until the solution *just* becomes cloudy. Clarify the solution by adding a few drops of the alcohol, and then determine its pH using pH paper.

Sodium Bicarbonate

Procedure

Dissolve a small amount of the compound in a 5% aqueous sodium bicarbonate solution. Observe the solution carefully. If the compound is an acid, you may see bubbles of carbon dioxide form. In some cases with solids, the evolution of carbon dioxide may not be that obvious.

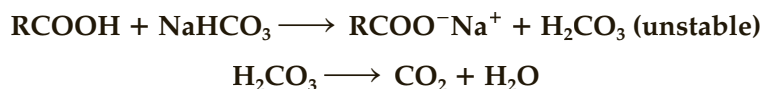

Silver Nitrate

Procedure

Acids may give a false silver nitrate test, as described in Experiment 57B.

Neutralization Equivalent (Optional)

Procedure

Accurately weigh (to three significant figures) approximately 0.2 g of the acid and place in a 125-mL Erlenmeyer flask. Dissolve the acid in about 50 mL of water or aqueous ethanol (the acid need not dissolve completely, because it will dissolve as it is titrated). Titrate the acid using a solution of sodium hydroxide of known molarity (about 0.1 M) and a phenolphthalein indicator.

Calculate the neutralization equivalent (NE) from the equation

$$\text{NE} = \frac{\text{mg acid}}{\text{molarity of NaOH} \times \text{mL of NaOH added}}$$

The NE is identical to the equivalent weight of the acid. If the acid has only one carboxyl group, the neutralization equivalent and the molecular weight of the acid are identical. If the acid has more than one carboxyl group, the neutralization equivalent equals the molecular weight of the acid divided by the number of carboxyl groups, that is, the equivalent weight. The NE can be used much like a derivative to identify a specific acid.

Some phenols are not sufficiently acidic to behave much like carboxylic acids. This is especially true of those substituted with electron-withdrawing groups at the *ortho* and *para* ring positions. These phenols, however, can usually be eliminated either by the ferric chloride test (see Experiment 57F) or by spectroscopy (phenols have no carbonyl group).

Spectroscopy

Infrared

C=O stretch is very strong and often broad in the region between 1725 cm⁻¹ and 1690 cm⁻¹.

O—H stretch has a very broad absorption in the region between 3300 cm⁻¹ and 2500 cm⁻¹; it usually overlaps the CH stretch region.


See Technique 25 for details.

Nuclear Magnetic Resonance




The acid proton of a —COOH group usually has resonance near 12.0 ppm. See Technique 26 for details. Carbon NMR is described in Technique 27.

Derivatives


Derivatives of acids are usually amides. They are prepared via the corresponding acid chloride:

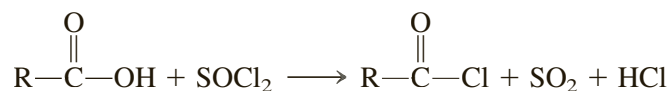

The most common derivatives are the amides, the anilides, and the *p*-toluidides.

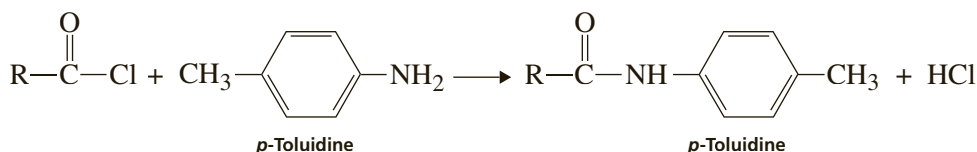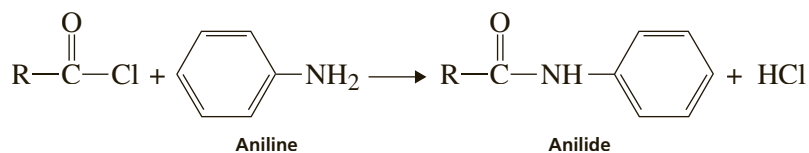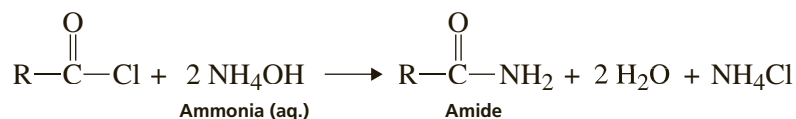

Procedures for the preparation of these derivatives are given in Appendix 2.

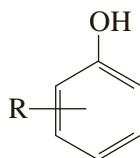57F EXPERIMENT 57F**Phenols**

Like carboxylic acids, phenols are acidic compounds. However, except for the nitrosubstituted phenols (discussed in the section covering solubilities), they are not as acidic as carboxylic acids. The $\text{p}K_a$ of a typical phenol is 10, whereas the $\text{p}K_a$ of a carboxylic acid is usually near 5. Hence, phenols are generally not soluble in the weakly basic sodium bicarbonate solution, but they dissolve in sodium hydroxide solution, which is more strongly basic.

Solubility Characteristics					Classification Tests
HCl	NaHCO ₃	NaOH	H ₂ SO ₄	Ether	Colored phenolate anion
(-)	(-)	(+)	(+)	(+)	Ferric chloride
Water: Most are insoluble, although phenol itself and the nitrophenols are soluble.					Ce(IV) Test Bromine/water

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous solutions in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container.

CLASSIFICATION TESTS

Sodium Hydroxide Solution

With phenols that have a high degree of conjugation possible in their conjugate base (phenolate ion), the anion is often colored. To observe the color, dissolve a small amount of the phenol in 10% aqueous sodium hydroxide solution. Some phenols do not give a color. Others have an insoluble anion and give a precipitate. The more acidic phenols, such as the nitrophenols, tend more toward colored anions.

Ferric Chloride

Procedure

Add about 50 mg of the unknown solid (2 mm or 3 mm off the end of a spatula) or 5 drops of the liquid unknown to 1 mL of water. Stir the mixture with a spatula so that as much as possible of the unknown dissolves in water. Add several drops of a 2.5% aqueous solution of ferric chloride to the mixture. Most water-soluble phenols produce an intense red, blue, purple, or green color. Some colors are transient, and it may be necessary to observe the solution carefully just as the solutions are mixed. The formation of a color is usually immediate, but the color may not last over any great period. Some phenols do not give a positive result in this test, so a negative test must not be taken as significant without other adequate evidence.

Test Compound

Try this test on phenol.

Discussion

The colors observed in this test result from the formation of a complex of the phenols with Fe(III) ion. Carbonyl compounds that have a high enol content also give a positive result in this test. The ferric chloride test works best with water-soluble phenols. A more reliable test, especially for water insoluble phenols, is the Ce(IV) test.

Cerium (IV) Test

Add 3 mL of 1,2-dimethoxyethane to 0.5 mL of Cerium(IV) reagent in a dry test tube. Gently shake the solution to thoroughly mix it, and then add 4 drops of a liquid compound to be tested. If you have a solid, you can directly add a few milligrams of the solid to the solution. Enough will dissolve to test if an —OH group is present. Gently shake the mixture, and look for an immediate color change from a yellow-orange solution to a red-orange or a deep red color indicating the presence of a phenol. The unsubstituted phenol, C₆H₅-OH, forms a dark-brown precipitate. Other phenols should yield a deep-red solution.

Test Compounds

Try this test on β -naphthol (2-naphthol).

Reagent

Prepare 2 M nitric acid solution by diluting 12.8 mL of concentrated nitric acid to 100 mL with water. Dissolve 8 g of ceric ammonium nitrate $[\text{Ce}(\text{NH}_4)_2(\text{NO}_3)_6]$ in 20 mL of the dilute nitric acid solution.

Discussion

The Ce(IV) test provides a more reliable way of detecting the presence of the hydroxyl group in water-insoluble phenols than the ferric chloride test. Since alcohols also give a color change with this reagent, you will first need to distinguish between alcohols and phenols by determining the solubility behavior of your compound. Phenols should be soluble in sodium hydroxide, whereas alcohols will not dissolve in aqueous sodium hydroxide.

Bromine Water**Procedure**


Prepare a 1% aqueous solution of the unknown, and then add a saturated solution of bromine in water to it, drop by drop while shaking, until the bromine color is no longer discharged. A positive test is indicated by the precipitation of a substitution product at the same time that the bromine color of the reagent is discharged.

Test Compound

Try this test on a 1% aqueous phenol solution.

Discussion

Aromatic compounds with ring-activating substituents give a positive test with bromine in water. The reaction is an aromatic substitution reaction that introduces bromine atoms into the aromatic ring at the positions *ortho* and *para* to the hydroxyl group. All available positions are usually substituted. The precipitate is the brominated phenol, which is generally insoluble because of its large molecular weight.

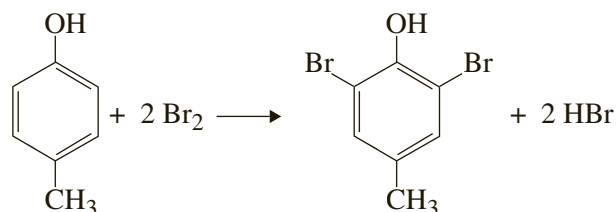

Other compounds that give a positive result with this test include aromatic compounds that have activating substituents other than hydroxyl. These compounds include anilines and alkoxyaromatics.

Spectroscopy**Infrared**

O—H stretch is observed near 3400 cm^{-1} .

C—O stretch is observed near 1200 cm^{-1} .

The typical aromatic ring absorptions between 1600 cm^{-1} and 1450 cm^{-1} are also found. Aromatic C—H is observed near 3100 cm^{-1} .

See Technique 25 for details.


Nuclear Magnetic Resonance

Aromatic protons are observed near 7 ppm. The hydroxyl proton has a resonance position that is concentration-dependent.

See Technique 26 for details. Carbon NMR is described in Technique 27.

Derivatives




Phenols form the same derivatives as alcohols (see Experiment 57H). They form urethanes on reaction with isocyanates. Phenylurethanes are used for alcohols, and the α -naphthylurethanes are more useful for phenols. Like alcohols, phenols yield 3,5-dinitrobenzoates.

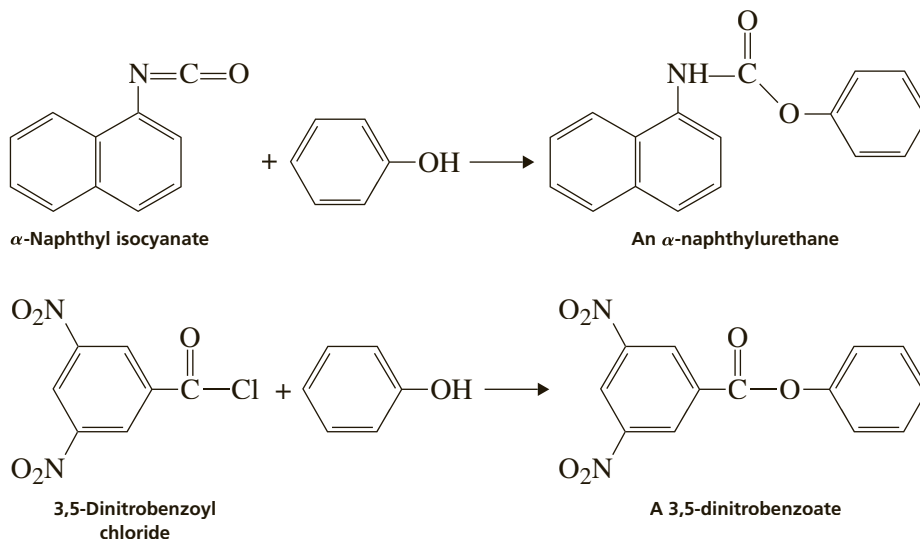

The bromine–water reagent yields solid bromo derivatives of phenols in several cases. These solid derivatives can be used to characterize an unknown phenol. Procedures for preparing these derivatives are given in Appendix 2.

57G EXPERIMENT 57G

Amines

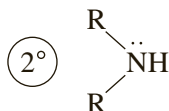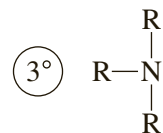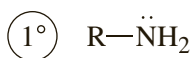

Amines are detected best by their solubility behavior and their basicity. They are the only basic compounds that will be issued for this experiment. Hence, once the

compound has been identified as an amine, the main problem that remains is to decide whether it is primary (1°), secondary (2°), or tertiary (3°). This can usually be decided either by the nitrous acid tests or by infrared spectroscopy.

Solubility Characteristics					Classification Tests
HCl	NaHCO ₃	NaOH	H ₂ SO ₄	Ether	pH of an aqueous solution
(+)	(-)	(-)	(+)	(+)	Hinsberg test
Water: < C ₆ (+)					Nitrous acid test
> C ₆ (-)					Acetyl chloride

SUGGESTED WASTE DISPOSAL

Residues from the nitrous acid test should be poured into a waste container containing 6 M hydrochloric acid. Dispose of all aqueous solutions in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container.

CLASSIFICATION TESTS

Nitrous Acid Test

Procedure

Dissolve 0.1 g of an amine in 2 mL of water to which 8 drops of concentrated sulfuric acid have been added. Use a large test tube. Often, a considerable amount of solid forms in the reaction of an amine with sulfuric acid. This solid is likely to be the amine sulfate salt. Add about 4 mL of water to help dissolve the salt. Any remaining solid will not interfere with the results of this test. Cool the solution to 5°C or less in an ice bath. Also cool 2 mL of 10% aqueous sodium nitrite in another test tube. In a third test tube, prepare a solution of 0.1 g β -naphthol in 2 mL of aqueous 10% sodium hydroxide, and place it in an ice bath to cool. Add the cold sodium nitrite solution, drop by drop while shaking, to the cooled solution of the amine. Look for bubbles of nitrogen gas. Be careful not to confuse the evolution of the *colorless* nitrogen gas with an evolution of *brown* nitrogen oxide gas. Substantial evolution of gas at 5°C or below indicates a primary aliphatic amine, RNH₂. The formation of a yellow oil or a yellow solid usually indicates a secondary amine, R₂NH. Either tertiary amines do not react, or they behave like secondary amines.

If little or no gas evolves at 5°C, take *half* the solution and warm it gently to about room temperature. Nitrogen gas bubbles at this elevated temperature indicate that the original compound was a **primary aromatic**, ArNH₂. Take the other half of the solution and, drop by drop, add the solution of β -naphthol in base. If a red dye precipitates, the unknown has been conclusively shown to be a primary aromatic amine, ArNH₂.

Test Compounds

Try this test with aniline, *N*-methylaniline, and butylamine.

CAUTION




The products of this reaction may include nitrosamines. Nitrosamines are suspected carcinogens. Avoid contact and dispose of all residues by pouring them into a waste container that contains 6 M hydrochloric acid.

Discussion

Before you make this test, it should definitely be proved by some other method that the unknown is an amine. Many other compounds react with nitrous acid (phenols, ketones, thiols, amides), and a positive result with one of these could lead to an incorrect interpretation.

The test is best used to distinguish *primary* aromatic and *primary* aliphatic amines from secondary and tertiary amines. It also differentiates aromatic and aliphatic primary amines. It cannot distinguish between secondary and tertiary amines. You will need to use infrared spectroscopy to make the distinction between secondary and tertiary amines. Primary aliphatic amines lose nitrogen gas at low temperatures under the conditions of this test. Aromatic amines yield a more stable diazonium salt and do not lose nitrogen until the temperature is elevated. In addition, aromatic diazonium salts produce a red azo dye when β -naphthol is added. Secondary and tertiary amines produce yellow nitroso compounds, which may be soluble or may be oils or solids. Many nitroso compounds have been shown to be carcinogenic. Avoid contact and immediately dispose of all such solutions in an appropriate waste container.

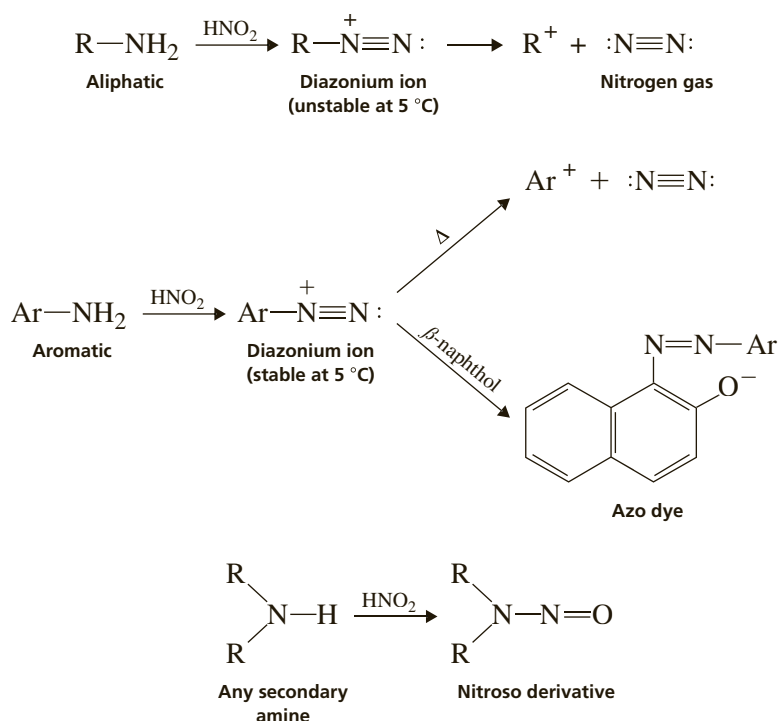

Hinsberg Test

A traditional method for classifying amines is the **Hinsberg test**. A discussion of this test can be found in the comprehensive textbooks listed prior to Section 57A. We have found that infrared spectroscopy is a more reliable method for distinguishing between primary, secondary, and tertiary amines.

pH of an Aqueous Solution

Procedure

If the compound is soluble in water, simply prepare an aqueous solution and check the pH with pH paper. If the compound is an amine, it will be basic and the solution will have a high pH. Compounds that are insoluble in water can be dissolved in ethanol–water or 1,2-dimethoxyethane–water.

Acetyl Chloride**Procedure**

Primary and secondary amines give a positive acetyl chloride test result (liberation of heat). This test is described for alcohols in Experiment 57H. Cautiously add dropwise the acetyl chloride to the liquid amine. This reaction can be very exothermic and violent! When the test mixture is diluted with water, primary and secondary amines often give a solid acetamide derivative; tertiary amines do not.

Test Compounds

Try this test with aniline and butylamine.

Spectroscopy**Infrared**

N—H stretch. Both aliphatic and aromatic primary amines show two absorptions (doublet due to symmetric and asymmetric stretches) in the region 3500–3300 cm^{-1} . Secondary amines show a single absorption in this region. Tertiary amines have no N—H bonds.

N—H bend. Primary amines have a strong absorption at 1640–1560 cm^{-1} . Secondary amines have an absorption at 1580–1490 cm^{-1} .

Aromatic amines show bands typical for the aromatic ring in the region 1600–1450 cm^{-1} .

Aromatic C—H is observed near 3100 cm^{-1} .

See Technique 25 for details.


Nuclear Magnetic Resonance

The resonance position of amino hydrogens is extremely variable. The resonance may also be very broad (quadrupole broadening). Aromatic amines give resonances near 7 ppm due to the aromatic ring hydrogens.


See Technique 26 for details. Carbon NMR is described in Technique 27.

Derivatives

The derivatives of amines that are most easily prepared are the acetamides and the benzamides. These derivatives work well for both primary and secondary amines but not for tertiary amines.

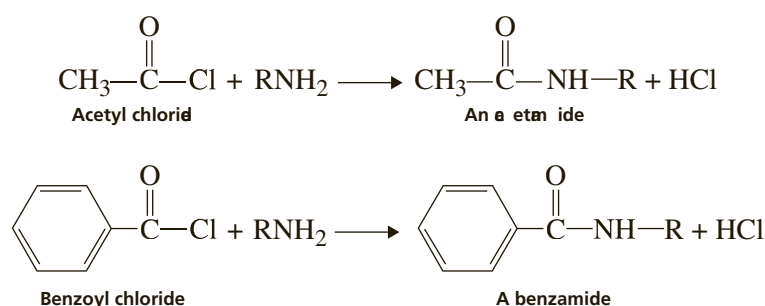

For tertiary amines, the methiodide salt is often useful.

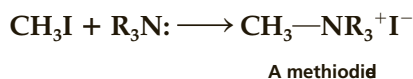


Procedures for preparing derivatives from amines can be found in Appendix 2.

57H EXPERIMENT 57H

Alcohols

Alcohols are neutral compounds. The only other classes of neutral compounds used in this experiment are the aldehydes, ketones, and esters. Alcohols and esters usually do not give a positive 2,4-dinitrophenylhydrazine test; aldehydes and ketones do. Esters do not react with Ce(IV) or acetyl chloride or with Lucas reagent, as alcohols do, and they are easily distinguished from alcohols on this basis. Primary and secondary alcohols are easily oxidized; esters and tertiary alcohols are not. A combination of the Lucas test and the chromic acid test will differentiate among primary, secondary, and tertiary alcohols.

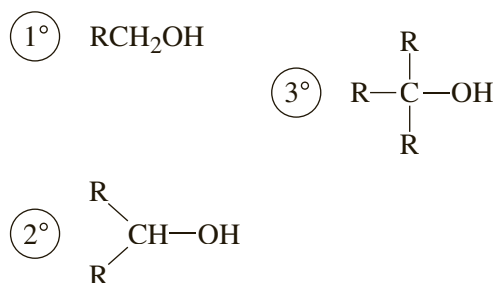

Solubility Characteristics					Classification Tests
HCl	NaHCO ₃	NaOH	H ₂ SO ₄	Ether	Cerium(IV) test
(-)	(-)	(-)	(+)	(+)	Acetyl chloride
Water: < C ₆ (+)					Lucas test
> C ₆ (-)					Chromic acid test
					Iodoform test

SUGGESTED WASTE DISPOSAL

Any solution containing chromium must be disposed of by placing it in a waste container specifically identified for the disposal of chromium wastes. Dispose of all other aqueous solutions in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container.

CLASSIFICATION TESTS

Cerium (IV) Test

Procedure for Water-Soluble or Partially Soluble Compounds

Add 3 mL of water to 0.5 mL of Cerium(IV) reagent in a test tube. Gently shake the solution to thoroughly mix it, and then add 4 drops of the compound to be tested. Gently shake the mixture and look for an *immediate color change* from a yellow-orange solution to a red-orange or deep red color indicating the presence of an —OH group in an alcohol or phenol. Phenol forms a dark-brown precipitate.

Test Compounds

Try this test on 1-butanol, 2-pentanol, 2-methyl-2-butanol, phenol, butanal, cyclohexanone, and ethyl acetate.

Procedure for Water-Insoluble Compounds

Add 3 mL of 1,2-dimethoxyethane to 0.5 mL of Cerium(IV) reagent in a dry test tube. Gently shake the solution to thoroughly mix it, and then add 4 drops of a liquid compound to be tested. If you have a solid, you can directly add a few milligrams of the solid to the solution. Enough will dissolve to test if an —OH group is present. Gently shake the mixture, and look for an *immediate color change* from a yellow-orange solution to a reddish-brown color indicating the presence of an alcohol or phenol.

Test Compounds

Try this test on 1-octanol, β -naphthol (2-naphthol), and benzoic acid.

Reagent

Prepare 2 M nitric acid solution by diluting 12.8 mL of concentrated nitric acid with 100 mL of water. Dissolve 8 g of ceric ammonium nitrate $[\text{Ce}(\text{NH}_4)_2(\text{NO}_3)_6]$ in 20 mL of the dilute nitric acid solution.

Discussion

Primary, secondary, and tertiary alcohols and phenols form 1:1 colored complexes with Ce(IV) and are an excellent way to detect hydroxyl groups. However, this is limited to compounds with no more than 10 carbon atoms. Unfortunately, the test cannot distinguish between primary, secondary, and tertiary alcohols. The Lucas test or chromium oxide test will have to be used for this purpose. Esters, ketones, carboxylic acids, and simple aldehydes do not change the color of the reagent and give a negative test with the Ce(IV) reagent. Thus, esters and other neutral compounds can be distinguished from alcohols by this test. Amines produce a flocculent white precipitate with the reagent. Cerium solutions can oxidize alcohols, but this usually occurs when the solution is heated or when the alcohol is in contact with the reagent for long periods.

Acetyl Chloride**Procedure**


Cautiously add about 5–10 drops of acetyl chloride, drop by drop, to about 0.25 mL of the liquid alcohol contained in a small test tube. Evolution of heat and hydrogen chloride gas indicates a positive reaction. Check for the evolution of HCl with a piece of wet blue litmus paper. Hydrogen chloride will turn the litmus paper red. Adding water will sometimes precipitate the acetate.

Test Compounds

Try this test with 1-butanol.

Discussion

Acid chlorides react with alcohols to form esters. Acetyl chloride forms acetate esters.

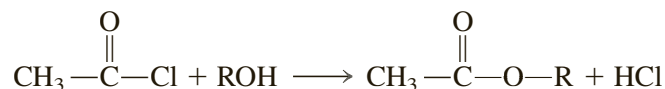

Usually, the reaction is exothermic, and the heat evolved is easily detected. Phenols react with acid chlorides somewhat as alcohols do. Hence, phenols should be

eliminated as possibilities before this test is attempted. Amines also react with acetyl chloride to evolve heat (see Experiment 57G). This test does not work well with solid alcohols.

Lucas Test

Procedure

Place 2 mL of Lucas reagent in a small test tube, and add 3–4 drops of the alcohol. Stopper the test tube and shake it vigorously. Tertiary (3°), benzylic, and allylic alcohols give an immediate cloudiness in the solution as the insoluble alkyl halide separates from the aqueous solution. After a short time, the immiscible alkyl halide may form a separate layer. Secondary (2°) alcohols produce a cloudiness after 2–5 minutes. Primary (1°) alcohols dissolve in the reagent to give a clear solution (no cloudiness). Some secondary alcohols may have to be heated slightly to encourage reaction with the reagent.

NOTE: This test works only for alcohols that are soluble in the reagent. This often means that alcohols with more than six carbon atoms cannot be tested.

Test Compounds

Try this test with 1-butanol (*n*-butyl alcohol), 2-butanol (*sec*-butyl alcohol), and 2-methyl-2-propanol (*t*-butyl alcohol).


Reagent

Cool 10 mL of concentrated hydrochloric acid in a beaker, using an ice bath. While still cooling and while stirring, dissolve 16 g of anhydrous zinc chloride in the acid.

This test depends on the appearance of an alkyl chloride as an insoluble second layer when an alcohol is treated with a mixture of hydrochloric acid and zinc chloride (Lucas reagent):



Primary alcohols do not react at room temperature; therefore, the alcohol is seen simply to dissolve. Secondary alcohols react slowly, whereas tertiary, benzylic, and allylic alcohols react instantly. These relative reactivities are explained on the same basis as the silver nitrate reaction, which is discussed in Experiment 57B. Primary carbocations are unstable and do not form under the conditions of this test; hence, no results are observed for primary alcohols.

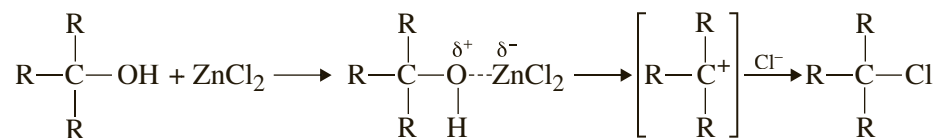

The Lucas test does not work well with solid alcohols or liquid alcohols containing six or more carbon atoms.

Chromic Acid Test: Alternative Test

CAUTION

Many chromium(VI) compounds are suspected carcinogens. If you would like to run this test, talk to your instructor first. The Lucas test will distinguish between 1°, 2°, and 3° alcohols, and you should do that test first. If you run the chromic acid test, be sure to wear gloves to avoid contact with this reagent.



Procedure

Dissolve 1 drop of a liquid or about 10 mg of a solid alcohol in 1 mL of *reagent-grade* acetone. Add 1 drop of the chromic acid reagent, and note the result that occurs within 2 seconds. A positive test for a primary or a secondary alcohol is the appearance of a blue-green color. Tertiary alcohols do not produce the test result within 2 seconds, and the solution remains orange. To make sure that the acetone solvent is pure and does not give a positive test result, add 1 drop of chromic acid to 1 mL of acetone that does not have an unknown dissolved in it. The orange color of the reagent should persist for *at least* 3 seconds. If it does not, a new bottle of acetone should be used.

Test Compounds


Try this test with 1-butanol (*n*-butyl alcohol), 2-butanol (*sec*-butyl alcohol), and 2-methyl-2-propanol (*t*-butyl alcohol).

Reagent


Dissolve 20 g of chromium trioxide (CrO₃) in 60 mL of cold water in a beaker. Add a magnetic stir bar to the solution. With stirring, slowly and carefully add 20 mL of concentrated sulfuric acid to the solution. This reagent should be prepared fresh each term.

Discussion

This test is based on the reduction of chromium(VI), which is orange, to chromium(III), which is green, when an alcohol is oxidized by the reagent. A change in color of the reagent from orange to green represents a positive test. Primary alcohols are oxidized by the reagent to carboxylic acids; secondary alcohols are oxidized to ketones.

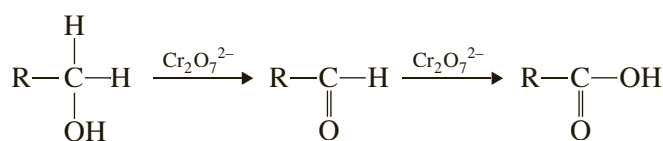

Primary alcohols

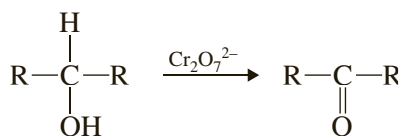

Secondary alcohols

Although primary alcohols are first oxidized to aldehydes, the aldehydes are further oxidized to carboxylic acids. The ability of chromic acid to oxidize aldehydes but not ketones is taken advantage of in a test that uses chromic acid to distinguish between aldehydes and ketones (see Experiment 57D). Secondary alcohols are oxidized to ketones, but no further. Tertiary alcohols are not oxidized at all by the reagent; hence, this test can be used to distinguish primary and secondary alcohols from tertiary alcohols. Unlike the Lucas test, this test can be used with all alcohols regardless of molecular weight and solubility.

Iodoform Test

Alcohols with the hydroxyl group at the 2-position of the chain give a positive iodoform test. See the discussion in Experiment 57D.

Spectroscopy*Infrared*

O—H stretch. A medium to strong, and usually broad, absorption comes in the region $3600\text{--}3200\text{ cm}^{-1}$. In dilute solutions or with little hydrogen bonding, there is a sharp absorption near 3600 cm^{-1} . In more concentrated solutions, or with considerable hydrogen bonding, there is a broad absorption near 3400 cm^{-1} . Sometimes both bands appear.

C—O stretch. There is a strong absorption in the region $1200\text{--}1500\text{ cm}^{-1}$. Primary alcohols absorb nearer 1050 cm^{-1} ; tertiary alcohols and phenols absorb nearer 1200 cm^{-1} . Secondary alcohols absorb in the middle of this range.

See Technique 25 for details.


Nuclear Magnetic Resonance

The hydroxyl resonance is extremely concentration-dependent, but it is usually found between 1 ppm and 5 ppm. Under normal conditions, the hydroxyl proton does not couple with protons on adjacent carbon atoms.


See Technique 26 for details. Carbon NMR is described in Technique 27.

Derivatives

The most common derivatives for alcohols are the 3,5-dinitrobenzoate esters and the phenylurethanes. Occasionally, the α -naphthylurethanes (Experiment 57F) are also prepared, but these latter derivatives are more often used for phenols.

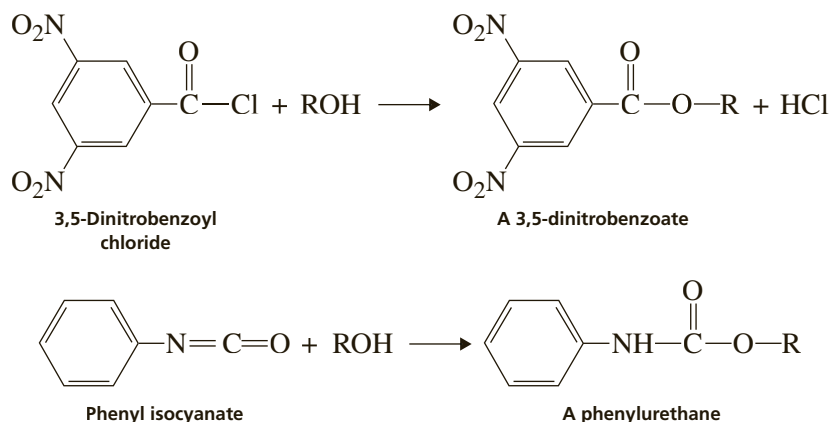

Procedures for preparing these derivatives are given in Appendix 2.

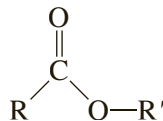57I EXPERIMENT 57I*Esters*

Esters are formally considered “derivatives” of the corresponding carboxylic acid. They are frequently synthesized from the carboxylic acid and the appropriate alcohol:



Thus, esters are sometimes referred to as though they were composed of an acid part and an alcohol part.

Although esters, like aldehydes and ketones, are neutral compounds that have a carbonyl group, they do not usually give a positive 2,4-dinitrophenylhydrazine test result. The two most common tests for identifying esters are the basic hydrolysis and ferric hydroxamate tests.

Solubility Characteristics					Classification Tests
HCl	NaHCO ₃	NaOH	H ₂ SO ₄	Ether	Ferric hydroxamate test
(-)	(-)	(-)	(+)	(+)	Basic hydrolysis
Water: < C ₄ (+)					
> C ₅ (-)					

SUGGESTED WASTE DISPOSAL

Solutions containing hydroxylamine or derivatives formed from it should be placed in a beaker containing 6 M hydrochloric acid. Dispose of any other aqueous solutions in the container designated for aqueous waste. Any remaining organic compounds must be disposed of in the appropriate organic waste container.

CLASSIFICATION TESTS

Ferric Hydroxamate Test

Procedure

Before starting, you must determine whether the compound to be tested already has enough enolic character in acid solution to give a positive ferric chloride test. Dissolve 1 or 2 drops of the unknown liquid or a few crystals of the unknown solid in 1 mL of 95% ethanol, and add 1 mL of 1 M hydrochloric acid. Add 1 or 2 drops of 5% ferric chloride solution. If a burgundy, magenta, or reddish-brown color appears, the ferric hydroxamate test cannot be used. It contains enolic character (see Experiment 57F).


If the compound did not show enolic character, continue as follows. Dissolve 5 or 6 drops of a liquid ester, or about 40 mg of a solid ester, in a mixture of 1 mL of 0.5 M hydroxylamine hydrochloride (dissolved in 95% ethanol) and 0.4 mL of 6 M sodium hydroxide. Heat the mixture until it boils for a few minutes. Cool the solution and then add 2 mL of 1 M hydrochloric acid. If the solution becomes cloudy, add 2 mL of 95% ethanol to clarify it. Add a drop of 5% ferric chloride solution, and note whether a color is produced. If the color fades, continue to add ferric chloride until the color persists. A positive test should give a deep burgundy, magenta, or reddish-brown color.

Test Compound


Try this test with ethyl butanoate.

Discussion

On being heated with hydroxylamine, esters are converted to the corresponding hydroxamic acids.

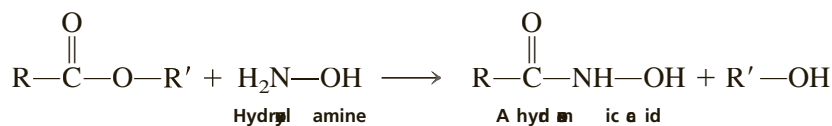

The hydroxamic acids form strong, colored complexes with ferric ion.

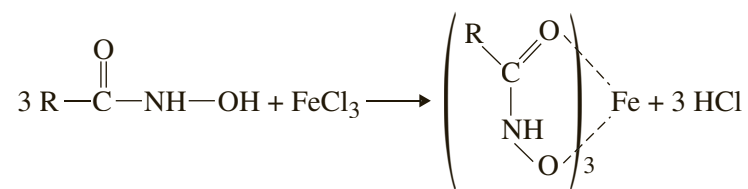

Basic Hydrolysis (Optional)


Procedure

Place 0.7 g of the ester in a 10-mL round-bottom flask with 7 mL of 25% aqueous sodium hydroxide. Add a boiling stone and attach a water condenser. Use a small amount of stopcock grease to lubricate the ground-glass joint. Boil the mixture for about 30 minutes. Stop the heating, and observe the solution to determine whether the oily ester layer has disappeared or whether the odor of the ester (usually pleasant) has disappeared. Low-boiling esters (below 110°C) usually dissolve within 30 minutes if the alcohol part has a low molecular weight. If the ester has not dissolved, reheat the mixture to reflux for 1–2 hours. After that time, the oily ester layer should have disappeared, along with the characteristic odor. Esters with boiling points up to 200°C should hydrolyze during this time. Compounds remaining after this extended period of heating are either unreactive esters or are not esters at all.

For esters derived from solid acids, the acid part can, if desired, be recovered after hydrolysis. Extract the basic solution with ether to remove any unreacted ester (even if it appears to be gone), acidify the basic solution with hydrochloric acid, and extract the acidic phase with ether to remove the acid. Dry the ether layer over anhydrous sodium sulfate, decant, and evaporate the solvent to obtain the parent acid from the original ester. The melting point of the parent acid can provide valuable information in the identification process.

Discussion

This procedure converts the ester to its separate acid and alcohol parts. The ester dissolves because the alcohol part (if small) is usually soluble in the aqueous medium, as is the sodium salt of the acid. Acidification produces the parent acid.

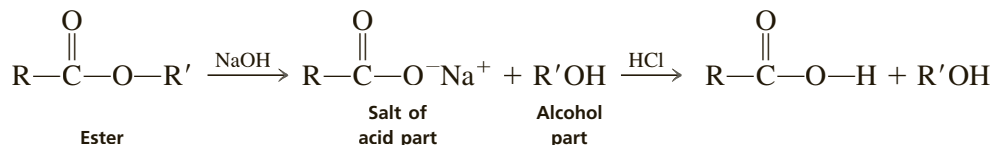

All derivatives of carboxylic acids are converted to the parent acid on basic hydrolysis. Thus, amides, which are not covered in this experiment, would also dissolve in this test, liberating the free amine and the sodium salt of the carboxylic acid.

Spectroscopy

Infrared

The ester-carbonyl group (C=O) peak usually indicates a strong absorption, as does the absorption of the carbonyl-oxygen link (C—O) to the alcohol part. C=O stretch at approximately 1735 cm⁻¹ is normal.¹ C—O stretch usually gives two or more absorptions, one stronger than the others, in the region 1280–1051 cm⁻¹.

See Technique 25 for details.

¹Conjugation with the carbonyl group moves the carbonyl absorption to lower frequencies. Conjugation with the alcohol oxygen raises the carbonyl absorption to higher frequencies. Ring strain (lactones) moves the carbonyl absorption to higher frequencies.

Nuclear Magnetic Resonance

Hydrogens that are alpha to an ester carbonyl group have resonance in the region 2–3 ppm. Hydrogens alpha to the alcohol oxygen of an ester have resonance in the region 2–5 ppm.


See Technique 26 for details. Carbon NMR is described in Technique 27.

Derivatives


Esters present a double problem when trying to prepare derivatives. To characterize an ester completely, you need to prepare derivatives of *both* the acid part and the alcohol part.

Acid Part


The most common derivative of the acid is the *N*-benzylamide derivative.

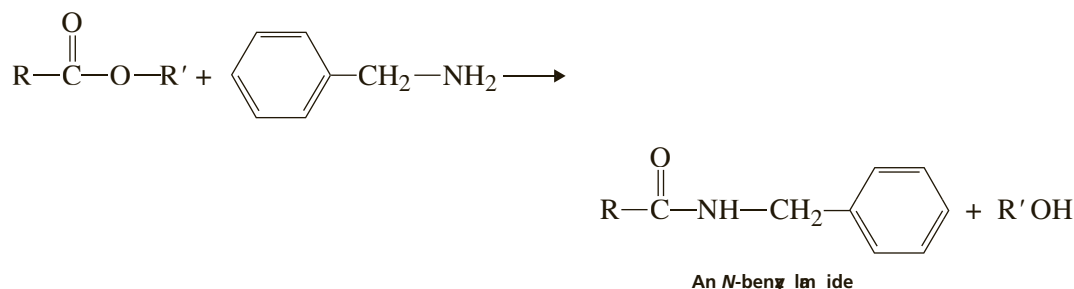

The reaction does not proceed well unless R⁺ is methyl or ethyl. For alcohol portions that are larger, the ester must be transesterified to a methyl or an ethyl ester before preparing the derivative.

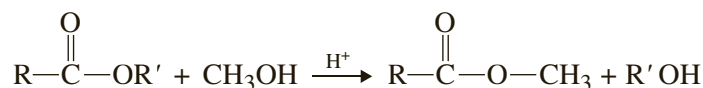


Hydrazine also reacts well with methyl and ethyl esters to give acid hydrazides.

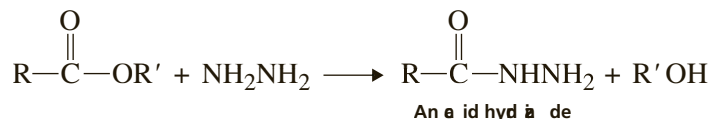

Alcohol Part

The best derivative of the alcohol part of an ester is the 3,5-dinitrobenzoate ester, which is prepared by an acyl interchange reaction.

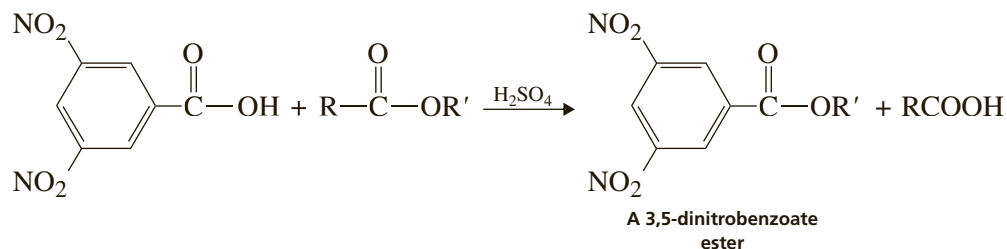

Most esters are composed of very simple acid and alkyl portions. For this reason, spectroscopy is usually a better method of identification than is the preparation of derivatives. Not only is it necessary to prepare two derivatives with an ester, but all esters with the same acid portion, or all those with the same alcohol portion, give identical derivatives of those portions.

Project-Based Experiments

Preparation of a C-4 or C-5 Acetate Ester

Esterification

Separatory funnel

Conventional distillation

In this experiment, we prepare an ester from acetic acid and a C-4 or a C-5 alcohol. This experiment is a conventional-scale preparation, but it is similar to the microscale preparation of isopentyl acetate, which is described in Experiment 14. However, for the experiment, either your instructor will assign, or you will pick, one of the following C-4 or C-5 alcohols to react with acetic acid:

1-Butanol (<i>n</i> -butyl alcohol)	1-Pentanol (<i>n</i> -pentyl alcohol)
2-Butanol (<i>sec</i> -butyl alcohol)	2-Pentanol
2-Methyl-1-propanol (isobutyl alcohol)	3-Pentanol
Cyclopentanol	3-Methyl-1-butanol (isopentyl alcohol)

If an NMR spectrometer is available, your instructor may wish to give you one of these alcohols as an unknown, leaving it to you to determine which alcohol was issued. For this purpose, you could use the infrared and NMR spectra, as well as the boiling points of the alcohol and its ester.

REQUIRED READING

Review: Essay Esters—Flavors and Fragrances
 Experiment 14
 Techniques 12, 13, and 14

SPECIAL INSTRUCTIONS

Be careful when dispensing sulfuric and glacial acetic acids. They are corrosive and will attack your skin if you make contact with them. If you get one of these acids on your skin, wash the affected area with large amounts of running water for 10–15 minutes.

If you select 2-butanol, reduce the amount of concentrated sulfuric acid to 0.5 mL. Also reduce the heating time to 60 minutes or less. Secondary alcohols have a tendency to give a significant percentage of elimination in strongly acidic solutions. Some of the alcohols may undergo elimination, leading to the formation of some low-boiling material (alkenes). In addition, cyclopentanol forms some dicyclopentyl ether, a solid.

SUGGESTED WASTE DISPOSAL

Any aqueous solutions should be placed in the container designated for dilute aqueous waste. Place any excess ester in the nonhalogenated organic waste container.

NOTES TO THE INSTRUCTOR

The sulfuric acid used as a catalyst in this reaction may be replaced with Dowex 50WX8-100 cationic exchange resin (sulfonate groups).

The purity of the esters can be determined by gas chromatography. It is recommended that a gas chromatogram of each of the starting alcohols be performed prior to determining the gas chromatogram of the esters. In this way, the peak corresponding to the parent alcohol can be identified by its retention time and the percentage of unreacted alcohol in the sample can be obtained. Approximate gas chromatography conditions for a GowMac Series 580 instrument with an 1/8-inch OV-1 column: 0.5 μL sample; flow rate, 27 mL/min; column temperature, 82°C; injector temperature, 170°C; detector temperature, 180°C; detector current, 200mA.

PROCEDURE

Apparatus

Assemble a reflux apparatus on top of your hot plate using a 20- or 25-mL round-bottom flask and a water-cooled condenser (refer to Figure 7.6A, but use a round-bottom flask instead of the conical vial). To control vapors, place a drying tube packed with calcium chloride on top of the condenser. Use a hot plate and the aluminum block with the larger set of holes for heating.

Reaction Mixture

Weigh (tare) an empty 10-mL graduated cylinder and record its weight. Place approximately 5.0 mL of your chosen alcohol in the graduated cylinder and reweigh it to determine the weight of alcohol. Disconnect the round-bottom flask from the reflux apparatus and transfer the alcohol into it. Do not clean or wash the graduated cylinder. Using the same graduated cylinder, measure approximately 7.0 mL of glacial acetic acid ($MW = 60.1$, $d = 1.06$ g/mL) and add it to the alcohol already in the flask. Using a calibrated Pasteur pipette, add 1 mL of concentrated sulfuric acid (0.5 mL if you have chosen 2-butanol), mixing *immediately* (swirl), to the reaction mixture contained in the flask. Add a corundum boiling stone or stirring bar and reconnect the flask. Do not use a calcium carbonate (marble) boiling stone, because it will dissolve in the acidic medium.

Reflux

Start water circulating in the condenser and bring the mixture to a boil. Continue heating under reflux for 60–75 minutes. Be sure to stir the mixture if you are using a stirring bar instead of a boiling stone. Then disconnect or remove the heating source and let the mixture cool to room temperature.

Extractions

Disassemble the apparatus and transfer the reaction mixture to a separatory funnel (60 or 125 mL) placed in a ring attached to a ring stand. Be sure the stopcock is

closed and, using a funnel, pour the mixture into the top of the separatory funnel. Also be careful to avoid transferring the boiling stone (or stirring bar), or you will need to remove it after the transfer. Add 10 mL of water, stopper the funnel, and mix the phases by careful shaking and venting (Section 12.7 and Figure 12.9). Allow the phases to separate and then uncap the funnel and drain the lower aqueous layer through the stopcock into a beaker or other suitable container. Next, extract the organic layer with 5 mL of 5% aqueous sodium bicarbonate just as you did previously with water. Extract the organic layer once again, this time with 5 mL of saturated aqueous sodium chloride.

Drying

Transfer the crude ester to a clean, dry, 25-mL Erlenmeyer flask and add approximately 1.0 g of anhydrous sodium sulfate. Cork the mixture and let it stand for 10–15 minutes while you prepare the apparatus for distillation. If the mixture does not appear dry (the drying agent clumps and does not “flow,” the solution is cloudy, or drops of water are obvious), transfer the ester to a new, clean, dry, 25-mL Erlenmeyer flask and add a new 0.5-g portion of anhydrous sodium sulfate to complete the drying.

Distillation

Assemble a distillation apparatus using your smallest round-bottom flask to distill from (Figure 14.10, but insert a water condenser as shown in Figure 14.11). Use a hot plate with an aluminum block to heat. Preweigh (tare) and use a 5-mL conical vial to collect the product. (It might be wise to have a second tared 5-mL conical vial handy in case you fill the first one.) Immerse the collection flask in a beaker of ice to ensure condensation and to reduce odors. If your alcohol is not an unknown, you can look up its boiling point in a handbook; otherwise, you can expect your ester to have a boiling point between 95 and 150°C. Continue distillation until only 1 or 2 drops of liquid remain in the distilling flask. Record the observed boiling point range in your notebook.

Yield Determination

Weigh the product and calculate the percentage yield of the ester. At the option of your instructor, determine the boiling point using one of the methods described in Technique 13, Section 13.2.

Spectroscopy

At your instructor’s option, obtain an infrared spectrum using salt plates (Technique 25, Section 25.2). Compare the spectrum with the one reproduced in Experiment 14. The spectrum of your ester should have similar features to the one shown. Interpret the spectrum and include it in your report to the instructor. You may also be required to determine and interpret the proton and carbon-13 NMR spectra (Technique 26, Sections 26.1 and Technique 27, Section 27.1). Submit your sample in a properly labeled vial with your report.

Gas Chromatography (Optional)

At your instructor’s option, perform a gas-chromatographic analysis of your ester. Either your instructor will provide a gas chromatogram of your starting alcohol or you will be asked to determine one at the same time that you do the analysis of your ester. Using both chromatograms, identify the alcohol and ester peaks and calculate the percentage of unreacted alcohol (if any) still remaining in your sample.

Is there any evidence of a product from a competing elimination reaction? Attach the chromatograms to your notebook or your final report, and be sure to include a discussion of the results in your report.

QUESTIONS

1. One method of favoring the formation of an ester is to add excess acetic acid. Suggest another method, involving the right-hand side of the equation that will favor the formation of the ester.
2. Why is the mixture extracted with sodium bicarbonate? Give an equation and explain its relevance.
3. Why are gas bubbles observed?
4. Using your alcohol, determine which starting material is the limiting reagent in this procedure. Which reagent is used in excess? How great is the molar excess (how many times greater)?
5. Outline a separation scheme for isolating your pure ester from the reaction mixture.
6. Interpret the principal absorption bands in the infrared spectrum of your ester or, if you did not determine the infrared spectrum of your ester, do this for the spectrum of isopentyl acetate (Experiment 14 and Technique 25 may be of some help).
7. Write a mechanism for the acid-catalyzed esterification that uses your alcohol and acetic acid. You may need to consult the chapter on carboxylic acids in your lecture textbook.
8. Tertiary alcohols do not react to produce an ester, like we would expect in this experiment. Instead, a different product is formed. When *t*-butyl alcohol (2-methyl-2-propanol) is used in this experiment, give the main organic product that would be formed.
9. Why is glacial acetic acid designated as “glacial”? (*Hint*: Consult a handbook of physical properties.)

Competing Nucleophiles in S_N1 and S_N2 Reactions: Investigations Using 2-Pentanol and 3-Pentanol

Nucleophilic substitution

Heating under reflux

Extraction

Gas chromatography

NMR spectroscopy

This experiment is based on the procedure outlined in Experiment 22. The purpose of this experiment is to examine the products formed when competing nucleophiles, equimolar concentrations of chloride ions and bromide ions, are allowed to react with either 2-pentanol or 3-pentanol. Based on the products formed in each reaction, students can advance a variety of hypotheses that account for the number and proportions of products formed.

Because the starting alcohols are secondary alcohols, one might expect that the substitution reactions will take place by a combination of S_N1 and S_N2 pathways. You will analyze the products of the three reactions in this experiment by a variety of techniques to determine the relative amounts of alkyl chloride and alkyl bromide formed in each reaction and to identify all of the products that are observed.

REQUIRED READING

Experiment 22	Nucleophilic Substitution Reactions: Competing Nucleophiles
Technique 7	Reaction Methods, Section 7.2, 7.4, 7.5, and 7.7
Technique 12	Extractions, Separations, and Drying Agents, Sections 12.5, 12.9, and 12.11
Technique 22	Gas Chromatography
Technique 26	Nuclear Magnetic Resonance Spectroscopy

Before beginning this experiment, review the appropriate chapters on nucleophilic substitution in your lecture textbook.

SPECIAL INSTRUCTIONS

Your instructor will also assign you either 2-pentanol or 3-pentanol. By sharing your results with other students, you will be able to collect data for both alcohols. To analyze the results of both experiments, your instructor will assign specific analysis procedures that the class will accomplish.

The solvent–nucleophile medium contains a high concentration of sulfuric acid. Sulfuric acid is corrosive; be careful when handling it. During the extractions, the longer your product remains in contact with water or aqueous sodium bicarbonate, the greater the risk that your product will decompose, leading to errors in your analytical results. Before coming to class, prepare so that you know exactly what you are supposed to do during the purification stage of the experiment.

SUGGESTED WASTE DISPOSAL

When you have completed the two experiments and all the analyses have been completed, discard any remaining alkyl halide mixture in the organic waste container marked for the disposal of halogenated substances. All aqueous solutions produced in this experiment should be disposed of in the container for aqueous waste.

NOTES TO THE INSTRUCTOR

The solvent–nucleophile medium must be prepared in advance for the entire class. Use the following procedure to prepare the medium. This procedure will provide enough solvent–nucleophile medium for about 10 students (assuming no spillage or other types of waste). Place 100 g of ice in a 500-mL Erlenmeyer flask and carefully add 76 mL concentrated sulfuric acid. Carefully weigh 19.0 g ammonium chloride and 35.0 g ammonium bromide into a beaker. Crush any lumps of the reagents to powder and then, using a powder funnel, transfer the halides to an Erlenmeyer flask. Carefully add the sulfuric acid mixture to the ammonium salts a little at a time. Swirl the mixture vigorously to dissolve the salts. It will probably be necessary to heat the mixture on a steam bath or a hot plate to achieve total solution. Keep a thermometer in the mixture and make sure that the temperature does not exceed 45°C. If necessary, you may add as much as 10 mL of water at this stage. Do not worry if a few small granules do not dissolve. When solution has been achieved, pour the solution into a container that can be kept warm until all students have taken their portions. The temperature of the mixture must be maintained at about 45°C to prevent precipitation of the salts. Be careful that the solution temperature does not exceed 45°C, however. Place a 20-mL calibrated pipette fitted with a pipette helper in the mixture. The pipette is always left in the mixture to keep it warm.

The gas chromatograph should be prepared as follows: Agilent (J & W) DB5 capillary column (30 m, 0.32 mm ID, 0.25 μm). Set injector temperature at 260°C. FID detector temperature is 280°C. The column oven conditions are as follows: start at 40°C (hold 2 min.), increase to 140°C at 20°C/min. (5 min.). The helium flow rate is 1.0 mL/min. The hydrogen FID gas makeup flow is 35 mL/min.

PROCEDURE

CAUTION



The solvent–nucleophile medium contains a high concentration of sulfuric acid. This liquid will cause severe burns if it touches your skin.

Apparatus

Assemble an apparatus for reflux using a 20-mL round-bottom flask, a reflux condenser, and a drying tube, as shown in the figure that accompanies Experiment 22. Loosely insert dry glass wool into the drying tube and then add water dropwise onto the glass wool until it is partially moistened. The moistened glass wool will trap the hydrogen chloride and hydrogen bromide gases produced during the reaction. As an alternative, you can use an external gas trap as described in Technique 7, Section 7.8, Part B. Do not place the round-bottom flask into the aluminum block until the reaction mixture has been added to the flask. Six Pasteur pipettes, two 3-mL conical vials with Teflon cap liners, and a 5-mL conical vial with a Teflon liner should also be assembled for use. All pipettes and vials should be clean and dry.

Preparation of Reagents

If a calibrated pipette fitted with a pipette helper is provided, you may adjust the pipette to 10 mL and deliver the solvent-nucleophile medium directly into your 20-mL round-bottom flask (temporarily placed in a beaker for stability). Alternatively, you may use a warm 10-mL graduated cylinder to obtain 10.0 mL of the solvent-nucleophile medium. The graduated cylinder must be warm in order to prevent precipitation of the salts. Heat it by running hot water over the outside of the cylinder or by putting it in the oven for a few minutes. Immediately pour the mixture into the round-bottom flask. With either method, a small portion of the salts in the flask may precipitate as the solution cools. Do not worry about this; the salts will redissolve during the reaction.

Reflux

Assemble the apparatus shown in the figure in Experiment 22. Using the following procedure, add 0.75 mL of either 2-pentanol or 3-pentanol, depending on which alcohol you were assigned, to the solvent-nucleophile mixture contained in the reflux apparatus. Dispense the alcohol from the automatic pipette or dispensing pump into a 10-mL beaker. Remove the drying tube and, with a 9-inch Pasteur pipette, dispense the alcohol directly into the round-bottom flask by inserting the Pasteur pipette into the opening of the condenser. Also add an inert boiling stone.¹ Replace the drying tube and start circulating the cooling water. Lower the reflux apparatus so that the round-bottom flask is in the aluminum block, as shown in the figure. Adjust the heat so that this mixture maintains a *gentle* boiling action. The aluminum block temperature should be about 140°C. Be careful to adjust the reflux ring, if one is visible, so that it remains in the lower fourth of the condenser. Violent boiling will cause loss of product. Heat the mixture for 45 minutes.

Purification

When the period of reflux has been completed, discontinue heating, lift the apparatus out of the aluminum block, and allow the reaction mixture to cool. Do not remove the condenser until the flask is cool. Be careful not to shake the hot solution as you lift it from the heating block or a violent boiling and bubbling action will result; this could allow material to be lost out of the top of the condenser. After the mixture has cooled for about 5 minutes, immerse the round-bottom flask (with condenser attached) in a beaker of cold tap water (no ice) and wait for this mixture to cool down to room temperature.

¹ Do not use calcium carbonate-based stones or Boileezers, because they will partially dissolve in the highly acidic reaction mixture.

There should be an organic layer present at the top of the reaction mixture. Add 0.75 mL of pentane to the mixture and *gently* swirl the flask. The purpose of the pentane is to increase the volume of the organic layer so that the following operations are easier to accomplish. Using a Pasteur pipette, transfer most (about 7 mL) of the bottom (aqueous) layer to another container. Be careful that the entire top organic layer remains in the boiling flask. Transfer the remaining aqueous layer and the organic layer to a 3-mL conical vial, taking care to leave behind any solids that may have precipitated. Allow the phases to separate and remove the bottom (aqueous) layer using a Pasteur pipette.

NOTE: For the following sequence of steps, be certain to be well prepared in advance. If you find that you are taking longer than 5 minutes to complete the entire extraction sequence, you probably have affected your results adversely!

Add 1.0 mL of water to the vial and *gently* shake this mixture. Allow the layers to separate and remove the aqueous layer, which is still on the bottom. Extract the organic layer with 1–2 mL of saturated sodium bicarbonate solution and remove the bottom aqueous layer.

Drying

Using a clean dry Pasteur pipette, transfer the remaining organic layer into a small test tube (10 × 75 mm) containing 3 to 4 microspatulafuls (using the V-grooved end) of anhydrous granular sodium sulfate. Stir the mixture with a microspatula, put a stopper on the tube, and set it aside for 10–15 minutes or until the solution is clear. If the mixture does not turn clear, add more anhydrous sodium sulfate. Transfer the halide solution with a clean, dry Pasteur pipette to a small, dry conical vial, taking care not to transfer any solid. The preferred method of storage is to use a Teflon stopper that has been secured *tightly* with a plastic cap. It is helpful to cover the stopper and cap with Parafilm (on the outside of the stopper and cap). Alternatively, you may use a screw-cap vial with a Teflon liner. *Be sure the cap is screwed on tightly.* Again, it is a good idea to cover the cap with Parafilm. Do not store the liquid in a container with a cork or a rubber stopper, because these will absorb the halides. If it is necessary to store the sample overnight, store it in a refrigerator. This sample can now be analyzed by as many of the methods as your instructor indicates.

ANALYSIS

The ratio of secondary pentyl chlorides and bromides must be determined. At your instructor's option, you may do this by gas chromatography or NMR spectroscopy or by both methods.

*Gas Chromatography*²

The instructor or a laboratory assistant may either make the sample injections or allow you to make them. In the latter case, your instructor will give you adequate

² *Note to the Instructor:* To obtain reasonable results for the gas-chromatographic analysis of the pentyl halides, it may be necessary to supply the students with response-factor corrections (Technique 22, Section 22.13). If pure samples of each product are available, check the assumption used here that the gas chromatograph responds equally to each substance. Response factors (relative sensitivities) are easily determined by injecting an equimolar mixture of products and comparing peak areas.

instruction beforehand. A reasonable sample size is 2.5 μL . Inject the sample into the gas chromatograph and record the gas chromatogram. The alkyl chlorides, because of their greater volatility, have a shorter retention time than the alkyl bromides.

Once the gas chromatogram has been obtained, determine the relative areas of the peaks (Technique 22, Section 22.12). If the gas chromatograph has an integrator, it will report the areas. Triangulation is the preferred method of determining areas if an integrator is not available. Record the percentages of all alkyl chloride and alkyl bromide products in the reaction mixture.

Nuclear Magnetic Resonance Spectroscopy

The instructor or a laboratory assistant will record the NMR spectrum of the reaction mixture.³ Submit a sample vial containing the mixture for this spectroscopic determination. The spectrum will also contain integration of the important peaks (Technique 26, "Nuclear Magnetic Resonance Spectroscopy"). Compare the integral of the *downfield* peaks of the alkyl halide multiplets. The relative heights of these integrals correspond to the relative amounts of each halide in the mixture.

REPORT

Record the percentages of all of the alkyl chloride and alkyl bromide products in the reaction mixture for each of the isomeric pentanol substrates. You will need to share your data with results obtained by someone in the class who used the other starting alcohol. Your laboratory report must include the percentages of each alkyl halide determined by each method used in this experiment for the two alcohols that were studied. On the basis of the products identified and their relative percentages, develop an argument for a mechanism that will account for all of the results obtained. All gas chromatograms and spectra should be attached to the report.

³ It is difficult to determine the ratio of pentyl chlorides and pentyl bromides using nuclear magnetic resonance. This method requires at least a 90-MHz instrument. At 300 MHz, the peaks are completely resolved.

Friedel–Crafts Acylation

Aromatic substitution

Directive groups

Vacuum distillation (optional)


Infrared spectroscopy

NMR spectroscopy (proton/carbon-13)

Structure proof

Gas chromatography (optional)

In this experiment, a Friedel–Crafts acylation of an aromatic compound is undertaken using acetyl chloride.

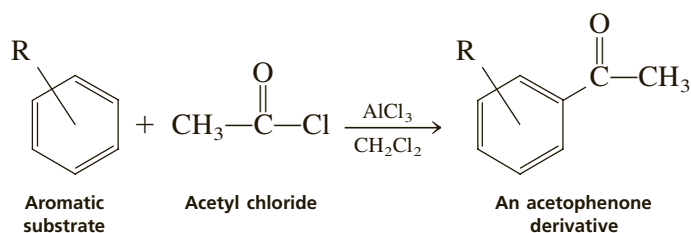

If benzene (R = H) were used as the substrate, the product would be a ketone, acetophenone. Instead of using benzene, however, you will perform the acylation on one of the following compounds:

Toluene

Ethylbenzene

o-Xylene (1,2-dimethylbenzene)

m-Xylene (1,3-dimethylbenzene)

p-Xylene (1,4-dimethylbenzene)

Pseudocumene (1,2,4-trimethylbenzene)

Anisole (methoxybenzene)

1,2-Dimethoxybenzene

1,3-Dimethoxybenzene

1,4-Dimethoxybenzene

Mesitylene (1,3,5-trimethylbenzene)

Hemellitrol

(1,2,3-trimethylbenzene, gives two products)

Except for the last one listed, each of these substrates will give a single product, a *substituted* acetophenone. You are to isolate this product and determine its structure by infrared and NMR spectroscopy. That is, you are to determine at which position of the original compound the new acetyl group becomes attached.

This experiment is much the same kind that professional chemists perform every day. A standard procedure, Friedel–Crafts acylation, is applied to a new compound for which the results are not known (at least not to you). A chemist who knows reaction theory well should be able to predict the result in each case. However, once the reaction is completed, the chemist must prove that the expected product has actually been obtained. If it has not, and sometimes surprises do occur, then the structure of the unexpected product must be determined.

To determine the position of substitution, several features of the product's spectra should be examined closely. These include the following:

INFRARED SPECTRUM

- **The C—H out-of-plane bending modes found between 900 and 690 cm^{-1} .** The C—H out-of-plane absorptions (see Technique 25, Figure 25.19A) often allow us to determine the type of ring substitution by their numbers, intensities, and positions.
- **The weak combination and overtone absorptions that occur between 2000 and 1667 cm^{-1} .** This set of combination bands (see Figure 25.19B) may not be as useful as the first set given here because the spectral sample must be very concentrated for them to be visible. But they are often weak. In addition, a broad carbonyl absorption may overlap and obscure this region, rendering it useless.

PROTON-NMR SPECTRUM

- **The integral ratio of the downfield peaks in the aromatic ring resonances found between 6 ppm and 8 ppm.** The acetyl group has a significant anisotropic effect, and those protons found *ortho* to this group on an aromatic ring usually have a greater chemical shift than the other ring protons (see Technique 26, Section 26.8 and Section 26.13).
- **A splitting analysis of the patterns found in the 6–8 ppm region of the NMR spectrum.** The coupling constants for protons in an aromatic ring differ according to their positional relations:

$$\textit{ortho} J = 6\text{--}10 \text{ Hz}$$

$$\textit{meta} J = 1\text{--}4 \text{ Hz}$$

$$\textit{para} J = 0\text{--}2 \text{ Hz}$$

If complex second-order splitting interaction does not occur, a simple splitting diagram will often suffice to determine the positions of substitution for the protons on the ring. For several of these products, however, such an analysis will be difficult. In other cases, one of the easily interpretable patterns like those described in Technique 26, Section 26.13 will be found.

CARBON-13 NMR SPECTRUM

- In *proton-decoupled* carbon-13 spectra, the number of resonances for the aromatic ring carbons (at about 120–130 ppm) will be of help in deciding the substitution patterns of the ring. Ring carbons that are equivalent by symmetry will give rise to a single peak, thereby causing the number of aromatic carbon peaks to fall below the maximum of six. A *p*-disubstituted ring, for instance, will show only four resonances. Carbons that bear a hydrogen usually will have a larger intensity than “quaternary” carbons. (See Technique 27, Section 27.6.)
- In *proton-coupled* carbon-13 spectra, the ring carbons that bear hydrogen atoms will be split into doublets, allowing them to be easily recognized.¹

¹ Note to the Instructor: For those not equipped to perform carbon-13 NMR spectroscopy, carbon-13 NMR spectra can be found reproduced in the Instructor's Manual.

As a final note, you should not eschew using the library. Technique 29 outlines how to find several important types of information. Once you think you know the identity of your compound, you might well try to find whether it has been reported previously in the literature and, if so, whether or not the reported data match your own findings. You may also wish to consult some spectroscopy books, such as Pavia, Lampman, Kriz, and Vyvyan, *Introduction to Spectroscopy*, or one of the other textbooks listed at the end of either Technique 25 or Technique 26, for additional help in interpreting your spectra.

REQUIRED READING

- Review:** Techniques 5, 6, 12, 25, 26, and 27
- | | |
|--------------|--|
| Technique 7 | Reaction Methods, Sections 7.5 and 7.8 |
| Technique 13 | Physical Constants of Liquids: The Boiling Point and Density |
- Optional:** Technique 16 Vacuum Distillation, Manometers, Sections 16.1, 16.2, and 16.8
- Optional:** Technique 22 Gas chromatography

Before you begin this experiment, you should review the chapters in your lecture textbook that deal with electrophilic aromatic substitution. Pay special attention to Friedel–Crafts acylation and to the explanations of directing groups. You should also review what you have learned about the infrared and NMR spectra of aromatic compounds.

SPECIAL INSTRUCTIONS

Both acetyl chloride and aluminum chloride are corrosive reagents. You should not allow them to come into contact with your skin, nor should you breathe them, because they generate HCl on hydrolysis. They may even react explosively on contact with water. When working with aluminum chloride, be especially careful to watch out for the powdered dust. Weighing and dispensing operations should be carried out in a hood. The workup procedure wherein excess aluminum chloride is decomposed with ice water should also be performed in the hood.

Your instructor will either assign you a compound or have you choose one for yourself from the list given at the beginning of this experiment. Although you will acetylate only one of these compounds, you should learn much more from this experiment by comparing your results with those of other students.

Notice that the details of vacuum distillation are left for you to figure out on your own. However, here are two hints. First, all of the products boil between 100°C and 150°C at 20-mm pressure. Second, if your chosen substrate is anisole, the product will be a *solid* with a low melting point and will solidify soon after vacuum distillation is completed. The solid can be distilled, but you should not run any cooling water through the condenser. It will also be worthwhile to preweigh the receiving flask, because it will be difficult to transfer the entire solidified product to another container to determine a yield.

SUGGESTED WASTE DISPOSAL

All aqueous solutions should be collected in a container specially marked for aqueous wastes. Place organic liquids in the container designated for nonhalogenated organic waste unless they contain methylene chloride. Waste materials that contain methylene chloride should be placed in the container designated for halogenated organic wastes. Note that your instructor may establish a different method of collecting wastes for this experiment.

NOTES TO THE INSTRUCTOR


It is suggested that you consider characterizing the Friedel-Crafts products by adding gas chromatography/mass spectrometry to the other spectroscopic techniques described in this experiment. Since most of the products show molecular ions, confirmation can be made of the molecular weight of the acylated product. The gas chromatography component will also help confirm that only a single acylated product was obtained. If the National Institute of Standards and Technology (NIST) database is available, confirmation of structure can be achieved.

You may want to consider omitting the vacuum distillation from this experiment. In almost all cases, a single product is formed, and the vacuum distillation does not materially improve the quality of the product. You may, however, observe unreacted starting material in the NMR spectrum and in the gas chromatographic analysis.

A four-step synthesis may be considered by linking together the Friedel-Crafts reaction with the synthesis of a chalcone (Experiment 65) and then preparing an epoxide (Experiment 66) from the chalcone and/or a cyclopropanated chalcone (Experiment 67). It is likely that the Friedel-Crafts reaction should produce enough acylated product for the reactions that follow. If you choose to link together the chalcone synthesis followed by epoxidation and cyclopropanation, it is suggested that you choose to prepare the acetyl derivatives of toluene, *p*-xylene, mesitylene, or anisole and use one of the recommended aldehydes shown in the following table to make the chalcone in Experiment 65.

Substrate	Aldehyde (Experiment 65)
toluene	4-methylbenzaldehyde
	4-nitrobenzaldehyde
	4-methoxybenzaldehyde
	piperonal
<i>p</i> -xylene	4-chlorobenzaldehyde
	4-fluorobenzaldehyde
	4-methoxybenzaldehyde
mesitylene	4-chlorobenzaldehyde
	4-methoxybenzaldehyde
anisole	4-chlorobenzaldehyde
	4-fluorobenzaldehyde
	4-methylbenzaldehyde
	piperonal

PROCEDURE

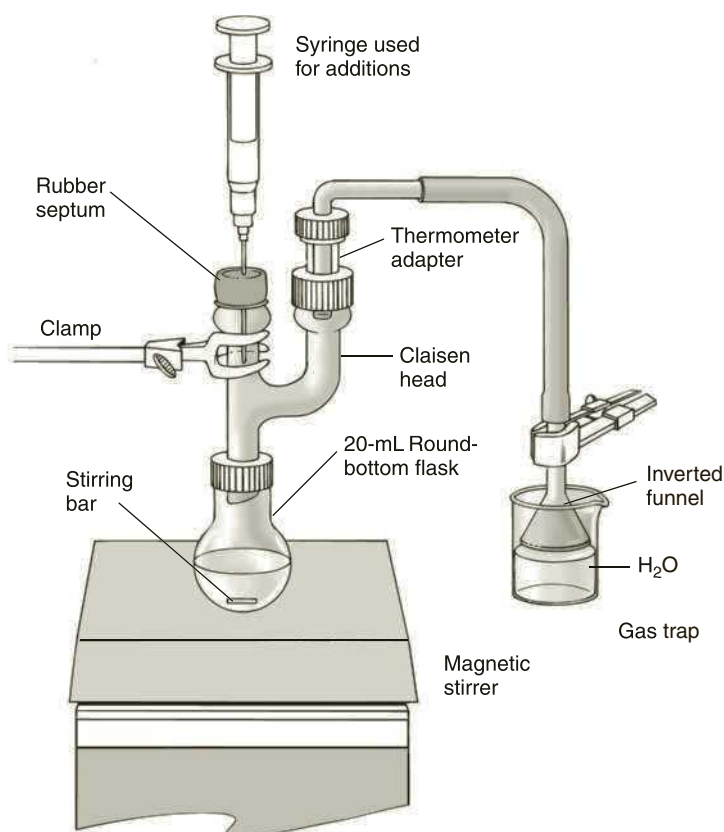

Apparatus for Friedel–Crafts reaction.

Assemble the reaction apparatus shown in the figure. It consists of a 20-mL round-bottom flask and a Claisen head with one opening fitted with a rubber septum and the other attached to an inverted-funnel trap for acidic gases. Secure the Claisen head and the gas-trap funnel with clamps. The funnel should be about 2 mm *above* the water. Remove the Claisen head and add 2 mL of methylene chloride, 0.8 g of AlCl_3 , and a magnetic stirring bar to the 20-mL round-bottom flask. Replace the Claisen head and begin stirring.

CAUTION



Both aluminum chloride and acetyl chloride are corrosive and noxious. Avoid contact and conduct all weighings in a hood. On contact with water, either compound may react violently.

Fill your 1-mL syringe (needle attached) with no less than 0.5 mL of fresh acetyl chloride. Insert the syringe through the rubber septum cap (see figure) and add the

acetyl chloride slowly over a 2-minute period. (Rapid addition of the acetyl chloride may cause foaming.) Using a graduated pipette and pipette pump, transfer exactly 0.5 mL of your chosen aromatic compound² to a preweighed 3-mL conical vial. Determine the weight of material delivered by weighing on a balance. Take up the aromatic compound with your syringe and slowly add it through the rubber septum over a 5-minute period. (This should not be done hastily, because the reaction is very exothermic; the mixture may boil up into the Claisen head.) When the aromatic compound has been added, rinse the vial with 1 mL of methylene chloride and, using the syringe, add this rinse to the reaction flask. Continue stirring for at least 30 minutes after the final addition has been made.

Isolation of Product

Remove the gas trap from the Claisen head and take the remaining apparatus, including the stirrer, to the hood. With your syringe, slowly add 4 mL of *ice-cold water* to the reaction mixture over a 5-minute period while stirring slowly.³ Next, add 4 mL of concentrated HCl with a Pasteur pipette and then stir the mixture vigorously with the magnetic stirrer until all the aluminum salts dissolve. At this point, discontinue stirring and allow the organic layer to separate. If the organic layer

² For substrates that are solids, weigh out 0.35g. After weighing, add a minimum amount of CH_2Cl_2 to dissolve the solid.

³ You may need to transfer all of the reaction mixture to a larger container if the flask is too full.

does not separate cleanly, add 0.5 mL of methylene chloride, stir again, and allow the organic layer to separate. You may have to add up to 1.5 mL of methylene chloride to induce the organic layer to separate cleanly.

Decant the entire mixture into a 15-mL centrifuge tube, leaving the stirring bar behind. Transfer the lower organic layer to a 5-mL conical vial with a filter-tip pipette. Avoid transferring any of the aqueous layer. If necessary, add a small amount of water and reparate the layers that have been transferred to the conical vial. If a significant amount of the original highly acidic aqueous layer is present, violent foaming will occur in the next step. Add about 1 mL of 5% sodium bicarbonate to the conical vial containing the organic layer. Cap the vial and shake it *gently*. Carefully vent the vial by loosening the cap and resealing it after a few moments. Repeat this mixing several times until the evolution of CO₂ is no longer apparent.

Transfer the organic layer to a dry 3-mL conical vial (5-mL if necessary) and add 3 to 4 microspatulafuls of anhydrous sodium sulfate (use the V-grooved end). Cap the vial and set it aside for 10–15 minutes while the liquid is dried. If the liquid appears cloudy, shake the vial several times during the drying period or add more sodium sulfate.

The product dissolved in the methylene chloride solution is likely to be highly colored at this point. Some of the color can be removed by employing column chromatography. Gently push a small amount of cotton into the constricted end of a Pasteur pipette. Add about 3 cm of alumina to the pipette. Remove the methylene chloride solution from the drying agent with a Pasteur pipette, and add it directly to the dry alumina contained in the chromatography column. Collect the eluent in a *preweighed and dry* conical vial. After collecting the liquid, add about 1 mL of fresh methylene chloride to the column and collect this eluent in the same conical vial. In a hood, place the vial in a hot water bath regulated to a temperature of about 40°C and direct a stream of air into the vial to evaporate the methylene chloride (Technique 7, Figure 7.17A). Do not rush this process. Allow the methylene chloride to be driven off completely. Monitor the evaporation by checking the volume markings on the side of the vial. When the volume is constant, the methylene chloride has been removed. If your instructor directs you to do a vacuum distillation, perform the optional procedure given in the next section; otherwise, skip to the “Boiling Point Determination and Spectroscopy” section. Weigh the conical vial after the methylene chloride has been removed and determine the percentage yield.

Vacuum Distillation (Optional)

If you are using a sand bath to heat, you should preheat it to about 165°C while assembling the apparatus. Assemble the apparatus *above* the sand bath; do not lower it into the sand bath until you are ready to distill. If you are using an aluminum block, preheating will not be necessary.

NOTE: Review Technique 16, Sections 16.1, 16.2, and 16.4, before proceeding.

Assemble an apparatus for vacuum distillation using an aspirator as shown in Technique 16, Figure 16.5. A manometer should be attached as shown in Technique 16, Figure 16.13. A piece of stainless steel sponge should be placed in the bottom portion of the neck of the Hickman still to protect the distilled product from any bumping action. Do not pack the stainless steel sponge too tightly. You may wish to preweigh the Hickman head (without the packing) to avoid having to transfer the product to determine the yield. This will be especially convenient if anisole was used as the substrate in the reaction. Using an *empty* conical vial, evacuate

the system and check for any leaks. When there are no significant leaks, add a spin vane to the 3-mL conical vial containing the product (methylene chloride removed). Attach the vial to the distillation apparatus and reestablish the vacuum.

If using a sand bath, lower the apparatus to begin the distillation and cover the sand bath with aluminum foil. If using an aluminum block, begin heating after lowering the apparatus. Adjust the spin vane to its maximum rate of spin. If boiling, bumping, or refluxing has not occurred after 3 minutes of heating, you may increase the heat. A sand bath or aluminum block temperature in the range of 165–200°C will be required, depending on your compound. Once the distillate begins to appear on the walls of the Hickman still, the distillation proceeds rapidly. When no liquid remains in the 3-mL vial or when liquid is no longer distilling, raise the apparatus immediately to discontinue the distillation. If you overheat the vial, it may crack. Turn the hot plate off. Allow the apparatus to cool to room temperature and then vent the system. Transfer the product to a preweighed storage container and determine its weight. (If you preweighed your Hickman still, remove the stainless steel sponge and transfer the Hickman still to a beaker for weighing.) Calculate the percentage yield.

Boiling-Point Determination and Spectroscopy

At the instructor's option, determine the boiling point of your product using the microboiling-point method (Technique 13, Section 13.2). Determine both the infrared and the NMR spectra (proton and carbon-13). The infrared spectra may be determined neat, using salt plates (Technique 25, Section 25.2), except for the product from anisole, which is a solid. For this product, one of the solution spectrum techniques (Technique 25, Section 25.6) should be used. The proton NMR spectra can be determined as described in Technique 26, Section 26.1. Deuteriochloroform is also an excellent solvent for all the carbon-13 samples as described in Technique 27, Section 27.1. Any residual methylene chloride appears at 5.3 ppm in the proton spectrum and at 54 ppm in the carbon spectrum.

REPORT

In the usual fashion, you should report the boiling point (or melting point) of your product, calculate the percentage yield, and construct a separation scheme diagram. You should also give the actual structure of your product. Include the infrared and NMR spectra and discuss carefully what you learned from each spectrum. If it did not help you determine the structure, explain why not. As many peaks as possible should be assigned on each spectrum and all important features explained, including the NMR splitting patterns, if possible.⁴ Discuss any literature you consulted and compare the reported results with your own.

Explain in terms of aromatic substitution theory why the substitution occurred at the position observed and why a single substitution product was obtained. Could you have predicted the result in advance?


⁴ Some starting material may be observed in the NMR spectrum. You should consult your instructor.

REFERENCE


Schatz, P. F. Friedel–Crafts Acylation. *J. Chem. Educ.* **1979**, 56 (Jul), 480.

QUESTIONS

1. The following are all relatively inexpensive aromatic compounds that could have been used as substrates in this reaction. Predict the product or products, if any, that would be obtained on acylation of each of them using acetyl chloride.

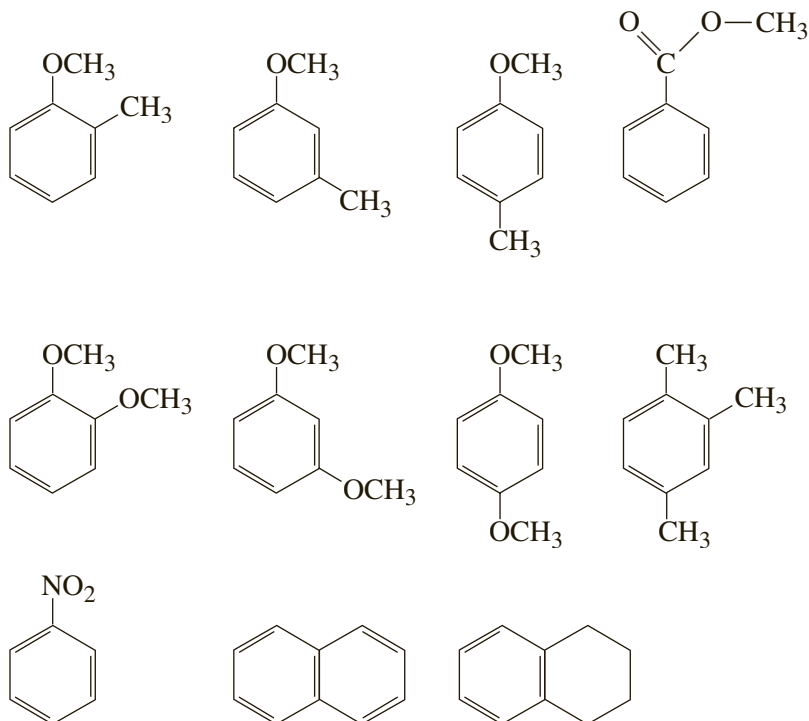

- Why are only monosubstitution products obtained in the acylation of the substrate compounds chosen for this experiment?
- Draw a full mechanism for the acylation of the compound you chose for this experiment. Include attention to any relevant directive effects.
- Why do none of the substrates given as choices for this experiment include any with meta-directing groups?
- Write equations for what happens when aluminum chloride is hydrolyzed in water. Do the same for acetyl chloride.
- Explain carefully, with a drawing, why the protons substituted ortho to an acetyl group normally have a greater chemical shift than the other protons on the ring.
- The compounds shown are possible acylation products from 1,2,4-trimethylbenzene (pseudocumene). Explain the only way you could distinguish these two products by NMR spectroscopy.

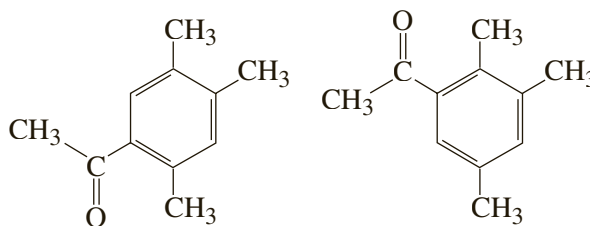

The Analysis of Antihistamine Drugs by Gas Chromatography–Mass Spectrometry

Gas chromatography–mass spectrometry

Critical thinking application

The use of **gas chromatography–mass spectrometry (GC–MS)** as an analytical technique is growing in importance. GC–MS is a powerful technique in which a gas chromatograph is coupled to a mass spectrometer that functions as the detector. If a sample is sufficiently volatile to be injected into a gas chromatograph, the mass spectrometer can detect each component in the sample and display its mass spectrum. The user can identify the substance by comparing its mass spectrum with the mass spectrum of a known substance. The instrument can also make this comparison internally by comparing the spectrum with spectra stored in its computer memory.

Antihistamines are a class of pharmaceutical agents commonly used to combat symptoms of allergies and colds. They reduce physiological effects of histamine production. Histamine, a protein, is normally released in the bloodstream as part of the body's reaction to intrusions by pollen, dust, molds, pet hair, and other **allergens** (substances that cause an allergic reaction). Even certain foods can cause an allergic response in some people. Excessive amounts of histamine can cause various disorders, including asthma, hay fever, sneezing, nasal secretions, skin irritations and swelling, hives, digestive disorders, and watering eyes. We take antihistamines to reduce these symptoms. Unfortunately, antihistamines also have some side effects, the most important of which is that they cause drowsiness. In fact, certain antihistamines are also sold as sleep aids.

In this experiment, you will prepare solutions of common over-the-counter antihistamine and cold-remedy tablets. The samples, once prepared, will be analyzed using a GC–MS instrument, and you will use the results to identify the significant antihistamine substances that comprise the original tablet.

REQUIRED READING

New: Technique 22 Section 22.12
Technique 28 Mass Spectrometry
Technique 29 Guide to the Chemical Literature

SPECIAL INSTRUCTIONS

This experiment requires the use of a GC–MS instrument. Before using this instrument, be certain to obtain instructions on its operation. As an option, your instructor may choose to perform the injections.

SUGGESTED WASTE DISPOSAL

Dispose of all solutions by discarding them in the container specified for nonhalogenated organic solvents. If your antihistamine contains either brompheniramine or chlorpheniramine, discard the solutions in the container specified for halogenated organic wastes.

PROCEDURE

Before beginning this experiment, you will need to rinse two 50-mL beakers, a syringe, and a snap-cap sample vial with HPLC-grade or spectrograde ethanol, and the glassware should be clean and dry before rinsing. Two rinsings of each item of glassware are recommended.

If your tablet has a colored coating, remove it by chipping it away from the tablet with a microspatula. Grind the tablet to a fine powder using a mortar and pestle. Weigh approximately 0.100 g of the powdered tablet into a prerinsed 50-mL beaker which has been prerinsed with ethanol. Add 10 mL of HPLC-grade ethanol to the beaker and let this solution stand, covered, for several minutes. Filter this solution by gravity through a fluted filter into a second prerinsed 50-mL beaker.

Draw the filtered solution into a prerinsed 5-mL syringe (without a needle), attach a 0.45- μm filter cartridge to the syringe, and expel the solution through the filter cartridge into a prerinsed sample vial. Repeat this process with a second syringe of solution. Cover the top of the sample vial with a square of aluminum foil and attach the snap-cap to the vial, over the top of the foil. Label the vial and store it in the refrigerator.

Analyze the sample by gas chromatography–mass spectrometry. Your instructor or laboratory assistant may either make the sample injections or allow you to make them. In the latter case, your instructor will give you adequate instruction beforehand. A reasonable sample size is 2 μL . Inject the sample into the gas chromatograph and obtain the printout of the total ion chromatogram, along with the mass spectrum of each component. You should also obtain the results of a library search for each component.

The library search will give you a list of components detected in your sample and the retention time and relative area for each component. The results will also list possible substances that the computer has tried to match against the mass spectrum of each component. This list—often called a “hit list”—will include the name of each possible compound, its Chemical Abstracts Registry number (CAS number), and a “quality” (“confidence”) measure, expressed as a percentage. The “quality” parameter estimates how closely the mass spectrum of the substance on the “hit list” fits the observed spectrum of that component in the gas chromatogram.

In your report, you should identify each significant component in the sample and provide its name and structural formula. You may have to use the CAS number as a key to look up the complete name and structure of the compound

(Technique 29, Section 29.11). You may need to search a computerized database to get the necessary information, or you may be able to find it in the *Aldrich Catalog Handbook of Fine Chemicals*, issued by the Aldrich Chemical Company. Current issues of this catalogue include listings of substances by CAS number. In your report, you should also report the relative percentage of the substance in the tablet extract. Finally, your instructor may also ask you to include the “quality” parameter from the “hit list.” If possible, determine which components have antihistamine activity and which ones are present for another purpose. The *Merck Index* may provide this information.

The Use of Organozinc Reagents in Synthesis: An Exercise in Synthesis and Structure Proof by Spectroscopy

Organometallic reactions

Green chemistry


Extractions

Use of a separatory funnel

Gas chromatography

Spectroscopy

A “green chemistry” alternative to the Grignard reaction was introduced in Experiment 36. This was the preparation of an organozinc reagent that was then allowed to react with a carbonyl compound. The reaction results in the formation of a new carbon–carbon bond.

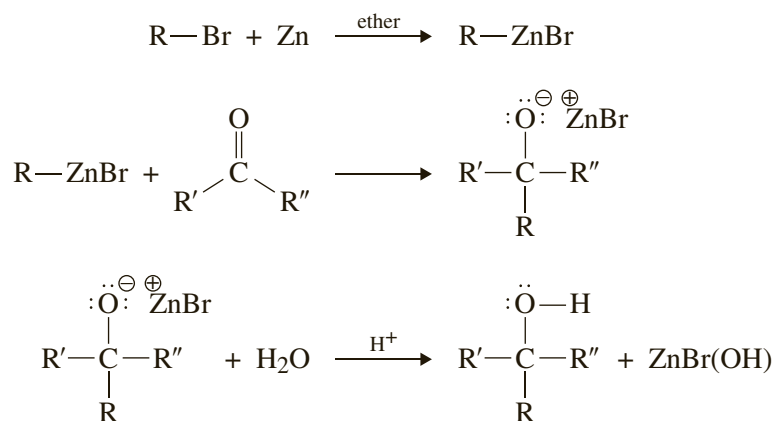

The theory that lies beneath this reaction was also presented.

The procedure that was outlined in Experiment 36 involved the preparation of a relatively simple alcohol. In this experiment, we will use the same general laboratory procedure, but we will use it to prepare substances that extend the synthetic methods presented in Experiment 36 and that will form products with interesting spectroscopic properties.

Each of the recommended starting carbonyl compounds has two different substituents attached to the carbonyl carbon. As a result, the products that are formed will contain a stereocenter (chiral carbon). The presence of the stereocenter will cause neighboring protons to become **diastereotopic** (see Technique 26, Section 26.16). Proton NMR spectra of these products will show the added complication resulting from unequal coupling to these diastereotopic protons. One of the purposes of this

experiment is to ask students to develop a complete analysis of the NMR spectra of the products formed in this reaction.

REQUIRED READING

<i>Review:</i>	Experiment 36	
	Technique 8	Section 8.3
	Technique 12	Sections 12.7, 12.8, 12.9, 12.11
	Technique 22	
	Technique 25	Sections 25.2, 25.4
	Technique 26	Section 26.1, 26.16
	Technique 27	Section 27.1

SPECIAL INSTRUCTIONS

This reaction involves the use of alkyl halides that are volatile and may also be **lachrymators**. Be certain to dispense these materials under the hood. Do not attempt to weigh these substances; determine the approximate volume of alkyl halide needed using the table of specific gravities provided in this experiment, and dispense the alkyl halides by volume using a calibrated pipette.

Students will be allowed to choose an alkyl halide and a carbonyl compound from the list provided.

SUGGESTED WASTE DISPOSAL

All aqueous solutions should be placed in a waste container designated for the disposal of aqueous wastes.

PROCEDURE

Activated Zinc

Carefully weigh 1.31 g (0.02 moles) of zinc powder and add it to a small (10-mL) Erlenmeyer flask or beaker. Add 1 mL of 5% aqueous hydrochloric acid and allow the mixture to stand for 1 to 2 minutes. There will be a noticeable evolution of hydrogen gas during this time. At the end of this period, pour the entire mixture into a Hirsch funnel and isolate the zinc by vacuum filtration. Rinse the zinc with 1 mL of water, followed by 1 mL of ethanol and 1 mL of diethyl ether. The zinc should be ready to use for the procedure, as described below.

Reaction of the Organozinc Reagent

Add 10 mL of saturated aqueous ammonium chloride solution to a 25-mL round-bottom flask. Add 1.31 g zinc powder (0.02 moles) and a stirring bar to the flask. Attach an air condenser to the flask and begin continuous stirring while adding the remaining reagents. Carefully weigh 0.01 moles of the carbonyl compound. Add the ketone or aldehyde and 1.6 mL of tetrahydrofuran to a test tube and add this solution dropwise to the zinc/ NH_4Cl solution. The rate of addition should be about 1 drop per second. Note that this addition can be made by dropping the solution


carefully down the opening in the air condenser; use a Pasteur pipette to add the solution. Allow the solution to stir for 10 to 15 minutes, giving time for the carbonyl compound to form a complex with the zinc. Add 0.02 moles (use the specific gravity to estimate the volume required) of the assigned alkyl halide to the stirring solution. *Be sure to dispense this reagent in the hood!* The rate of addition should be about 1 drop per second. Add the halide carefully by dropping it down the opening in the air condenser. Allow the reaction mixture to stir for 1 hour.

Set up a vacuum filtration with a Hirsch funnel. Decant the liquid from the reaction mixture through the Hirsch funnel. Rinse the round-bottom flask with approximately 1 mL of diethyl ether and pour the liquid into the Hirsch funnel. Using a second 1-mL portion of diethyl ether, rinse the solid that has collected in the Hirsch funnel. Discard the solid. Prepare a filter-tip pipette and transfer the liquid that was collected in the vacuum filtration into a separatory funnel. Use 1 mL of diethyl ether to rinse the inside of the filter flask and use the filter-tip pipette to transfer this liquid to the separatory funnel. Shake the separatory funnel gently to extract the organic material from the aqueous layer to the ether layer. Drain the lower (aqueous) layer into a 50-mL Erlenmeyer flask. *Do not discard this aqueous layer.* Collect the upper (organic) layer from the separatory funnel into a 25-mL Erlenmeyer flask (remember to collect the upper layer by pouring it from the top of the separatory funnel). Replace the aqueous layer in the separatory funnel and wash it with a 2-mL portion of ether. Separate the layers, save the aqueous layer in the same 50-mL Erlenmeyer flask as before, and combine the ether layer with the ether solution collected in the previous extraction. Repeat this extraction process of the aqueous phase one more time, using a fresh 2-mL portion of ether. Dry the combined ether extracts with 3–4 microspatulafuls of anhydrous sodium sulfate. Stopper the Erlenmeyer flask with a cork and allow it to stand for at least 15 minutes (or overnight).


Use a filter-tip pipette to transfer the dried liquid to a clean, preweighed Erlenmeyer flask. Use a small amount of ether to rinse the inside of the original flask and add this ether to the dried liquid. Evaporate the ether under a gentle stream of air. When the ether has evaporated completely, reweigh the flask to determine the yield of product. If it should be necessary to store your final product, use Parafilm to seal the container.

Prepare a sample of your final product for analysis by gas chromatography. Determine the infrared spectrum and both proton and ^{13}C NMR spectrum of your product. Use these spectra to determine the structure of your product. In your laboratory report, include an interpretation of each spectrum, identifying the principal absorption bands and demonstrating how the spectrum corresponds to the structure of your compound. Submit your sample in a labeled vial with your laboratory report.


Carbonyl starting materials

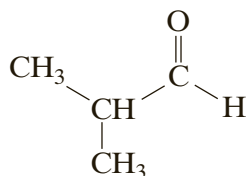

Isobutyraldehyde
MW = 72.1 g/mole

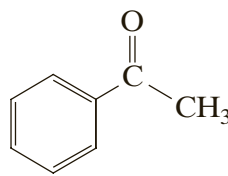


Acetophenone
MW = 120.2 g/mole

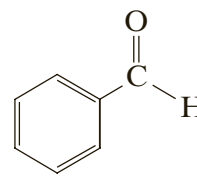

Benzaldehyde
MW = 106.1 g/mole


Alkyl halide starting materials

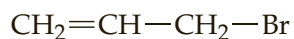

Allyl bromide
MW = 121.0 g/mole
Specific gravity = 1.398 g/mL



Crotyl chloride
MW = 90.6 g/mole
Specific gravity = 0.92 g/mL

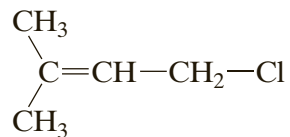

1-Chloro-3-methyl-2-butene
MW = 104.6 g/mole
Specific gravity = 0.98 g/mL

Synthesis of Naproxen by Palladium Catalysis

Preparation of an analgesic drug


Organometallic chemistry

Palladium-catalyzed reaction

Advanced laboratory experiment

Nonsteroidal anti-inflammatory drugs (NSAIDs) are among the most common drugs on the market. These drugs have analgesic (pain-reducing) and antipyretic (fever-reducing) effects. They also have anti-inflammatory effects. The common over-the-counter drugs in this group are ibuprofen, naproxen sodium, and aspirin. Naproxen sodium is most commonly sold as the over-the-counter drug, Aleve®. Originally, naproxen was synthesized by the Syntex Corporation, now owned by Roche Biosciences. They developed a process by which they could synthesize naproxen on a large scale.^{1,2} Aleve® is marketed by the Bayer HealthCare group.

Some NSAIDs are sold as racemic mixtures. One common example is ibuprofen, sold as Advil®, and marketed by Pfizer. Typically, only a single enantiomer is pharmacologically active. In the case of ibuprofen, it is the (*S*)-enantiomer. The inactive (*R*)-enantiomer in racemic ibuprofen is converted by an isomerase enzyme in the gut into the active enantiomer. Ibuprofen is now available as single enantiomer in certain preparations that supposedly has improved therapeutic value. Naproxen has always been marketed as the single enantiomer.

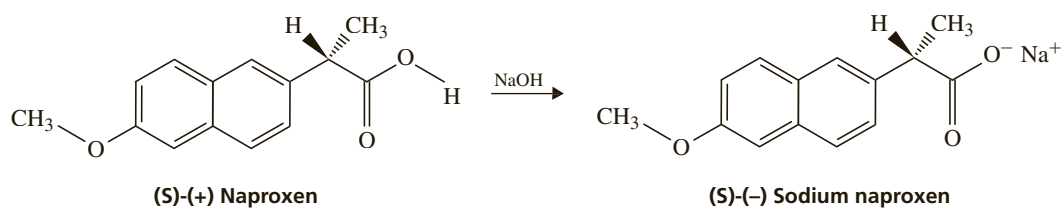

Notice in the structures shown above that the carboxylic acid, naproxen, has a plus sign of rotation while the sodium salt, sodium naproxen, has a minus sign of rotation when measured in a polarimeter. You should recall from your study of stereochemistry in your lecture course that a compound with (*S*) stereochemistry may have either a plus or a minus value in a polarimeter. Likewise, a compound with the (*R*) stereochemistry can have either a plus or minus value.

In this experiment, you will be preparing the carboxylic acid, naproxen, rather than the sodium salt, sodium naproxen. As is often the case with pharmaceuticals,

¹ <http://en.wikipedia.org/wiki/Naproxen>




² Harrington, P. J.; Lodewijk, E., Twenty years of Naproxen Technology, *Organic Process Research and Development* 1997, 1, 72–76.

the ionic compounds are found to be very soluble in an aqueous environment, thereby improving the transfer of the medicinal compound into the body. Naproxen, itself, would not be as soluble in an aqueous environment.

In this experiment we present a new process involving a palladium metal-catalyzed reaction sequence.^{3,4} While the process employed in this experiment is not competitive with the Roche/Syntex process, it provides an interesting use of an organometallic catalysis for instructional purposes.




The reaction involves the following sequence of reactions.

1. Butyllithium reacts with dicyclohexylamine to form lithium dicyclohexylamide (eq 1).
2. Lithium dicyclohexylamide, a very strong base, reacts with *tert*-butyl propionate to form the enolate (eq 2).
3. Pd⁰, the active form of the metal, is formed from the disproportionation of the palladium catalyst (eq 3).
4. Oxidative addition forms a complex of the aromatic halide with Pd⁰ to change the oxidation state of the palladium to Pd^{II} (eq 4).
5. The enolate replaces the halide in the complex. The oxidation state of Pd is still Pd^{II} (eq 5).
6. Reductive elimination forms *tert*-butyl naproxen, and regenerates the Pd⁰ active catalyst (eq 6).

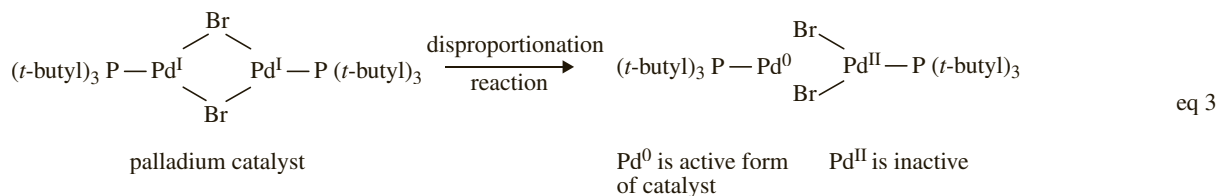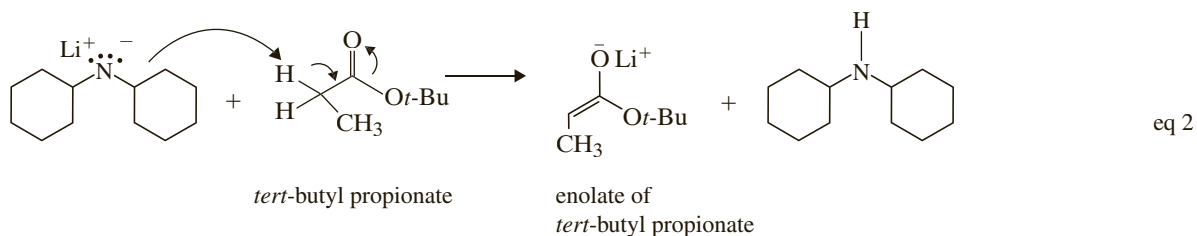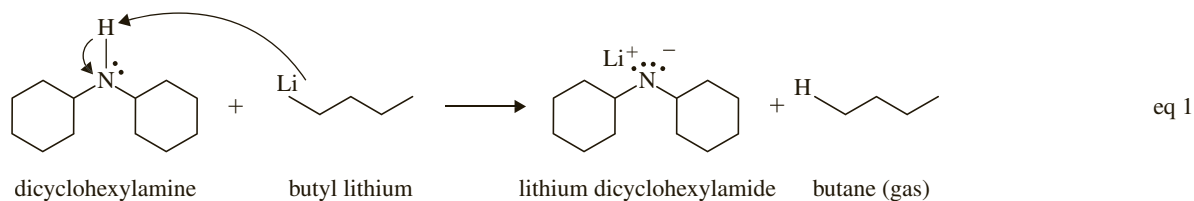


³ Hama, T.; Hartwig, J. F. *Org. Lett.* **2008**, *10*, 1545–1548.

⁴ Rumberger, S.; Visser, L.; Pittman, J.; Lampman, G. M., Department of Chemistry, Western Washington University

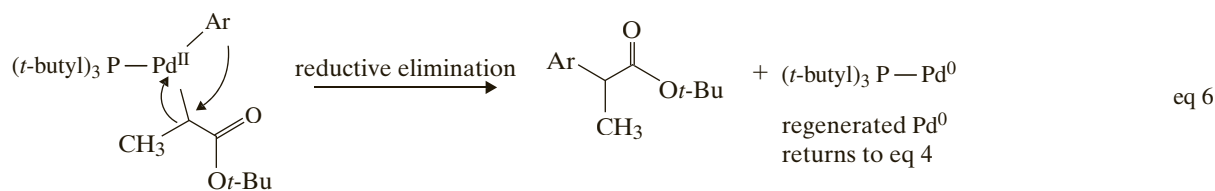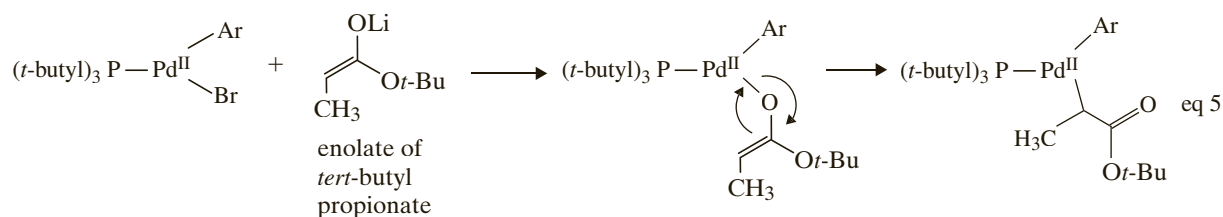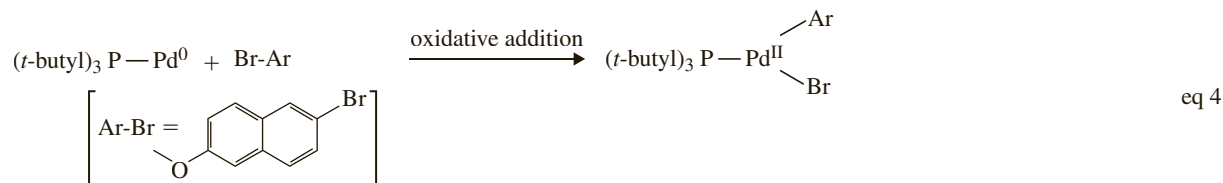

A simplified catalytic cycle is shown below. Notice that the formation of lithium dicyclohexylamide (eq 1) and the formation of the enolate (eq 2) are not included in the catalytic cycle.

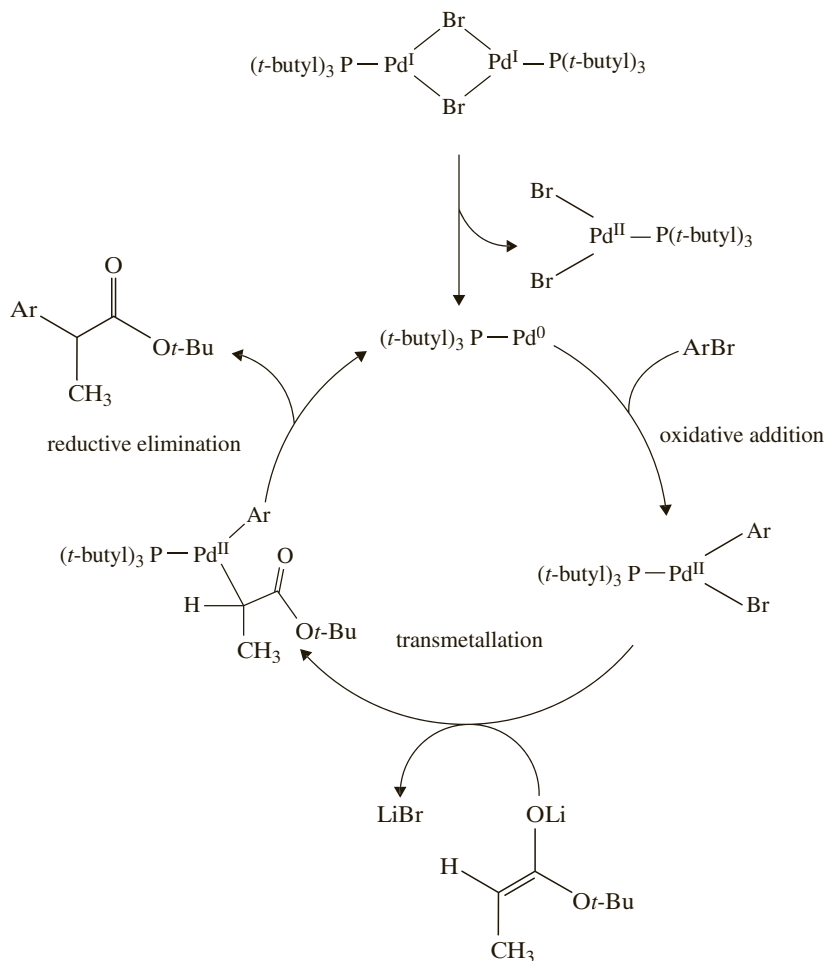

REQUIRED READING

Review: Techniques 5, 6, 7, 8, 10, 12, 21, 23, 25, and 26

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous solutions into the container designated for aqueous waste. Place the non-halogenated and halogenated organic waste into the appropriate waste containers.

NOTES TO THE INSTRUCTOR

This experiment is definitely more advanced than others in this book and should be carefully evaluated before including it in your organic laboratory program. You may want to consider it for a more advanced laboratory course.⁵ A pyrophoric compound, butyllithium, is used in this experiment. Pyrophoric compounds are known to catch fire when in the presence of air or moisture. Although 1.6 M butyllithium is among the safer organometallic reagents, *an instructor or trained assistant must be present at all times when using this reagent*. Since small amounts of butyllithium reagent are used, this helps to ensure safety. Butyllithium should be stored in a sealed system under an inert gas for reasons of safety. It is not nearly as dangerous to work with as the extremely pyrophoric *tert*-butyllithium. Nevertheless, it is good practice to use butyllithium with care under a dry inert atmosphere. <http://www.chemeurope.com/en/encyclopedia/N-Butyllithium.html>. The University of Minnesota guide to handling pyrophoric materials is very helpful. See pages 11–16 at the following website for a technique using a syringe to transfer butyllithium. http://www.dehs.umn.edu/PDFs/Pyrophoric_Chemicals_Guide.pdf. Syringes are required to maintain anhydrous conditions. The instructor's manual provides more information on handling butyllithium. The palladium catalyst should be stored in a desiccator or in a refrigerator. Before storing the catalyst, introduce some inert gas into the catalyst container, replace the cap, and completely cover the top with Parafilm.

The preparation of racemic naproxen can be completed in two laboratory periods and the experiment may be stopped here. It is suggested that you allow four laboratory periods in order to complete all of the parts of the experiment, including the resolution of the (*S*) enantiomers and the separation of the enantiomers using chiral HPLC. It is suggested that students work in pairs, to save time. An inert atmosphere, either nitrogen or argon, must be provided in order to maintain anhydrous and oxygen-free conditions. This experiment makes use of larger scale glassware including a 250-mL round-bottom flask with a F 19/22 joint and a 125-mL separatory funnel. A rotary evaporator should be made available in order to remove toluene solvent from the reaction mixture. Finally, although not required, it is helpful to have a chiral HPLC column available.

⁵ This experiment has been used successfully at Western Washington University for several years in the second quarter laboratory courses that meet 6 hours each week. The class includes chemistry majors, premedical, pre dental, and prepharmacy students.

PROCEDURE (WORK IN PAIRS)

Drying Glassware

Your instructor will place a large magnetic stir bar in a 250-mL round-bottom flask with a F 19/22 joint and dry the flask in an oven set to 125°C, for at least 3 hours, in advance of the beginning of the laboratory period. After the flask is dried in an oven, insert a #33 Suba Seal rubber septum and allow the flask to cool to room temperature.

Anhydrous Toluene

The instructor may provide 75 mL of anhydrous toluene dispensed from a commercial solvent purification system into a 250-mL round-bottom flask.⁶ An inert atmosphere remains in the flask when dispensed from a solvent purification system. Alternatively, anhydrous toluene may be purchased.⁷ In this case, you will need to purge the air from the flask containing toluene after it is dispensed. Once the toluene is dispensed it is essential to make sure that the dry conditions are rigorously maintained.

Maintaining Anhydrous Conditions

A gas manifold with multiple stopcocks and a single tank of a nitrogen or argon gas is an excellent way of running this reaction with a large laboratory class (see Technique 7, Figure 7.1 1B).⁸ This setup allows the contents of the flask to be maintained under an inert atmosphere. Attach syringe needles to each stopcock with a section of rubber tubing. Into each end of the rubber tubing, insert the barrel of a cut-off section of a plastic 1-mL syringe with an attached syringe needle. In advance, the instructor will turn on the inert gas source, and allow the gas to flow rapidly for 10 or 15 minutes to expel oxygen from the manifold and stopcocks and replace it with the inert gas. Open the stopcock on the manifold and insert the needle into the #33 Suba Seal rubber septum that is attached to the 250-mL round-bottom flask as shown in Technique 7, Figure 7.11B. Then, insert another syringe needle into the septum to provide an exit for the inert gas from the flask. Adjust the inert gas flow so that there are only a few bubbles emerging from the bubbler attached to the *exit port* shown in Figure 7.11B. With the stopcock opened, any residual oxygen that is present in the flask will be removed. This process maintains an inert gas atmosphere in the reaction flask.

Preparation of Lithium Dicyclohexylamide (eq 1)

Draw up 1.0 mL of dicyclohexylamine ($\text{C}_{12}\text{H}_{23}\text{N}$, MW 181 g/mol, 0.992 g, 0.0055 mol)⁹ into a 1-mL syringe with an attached needle and insert the needle into the rubber septum on the 250-mL flask. Add the amine from the syringe.¹⁰ Withdraw the empty syringe and needle after adding the amine. Start the stirrer.

⁶ Commercial solvent purification systems are often found in research laboratories in larger institutions. Since there are no graduations on the 250-mL round-bottom flask, take another 250-mL round-bottom flask, add 75 mL of water, and make a mark on the round-bottom flask to indicate the approximate location of the 75-mL mark. Then make a similar mark on the pre-dried 250-mL round-bottom flask, so that the student will be able to know how much toluene is to be added from the solvent purification system.

⁷ Aldrich Chemical Company, anhydrous toluene, #244511 sold in Sure/Seal containers to maintain anhydrous conditions.

⁸ It is advised to use a single gas manifold, rather than a double one. If a manifold is not available, an alternative method is described in the instructor's manual.

⁹ Aldrich Chemical Company, dicyclohexylamine 99%, #185841, used as supplied without purification.

¹⁰ It may be helpful to preweigh the syringe, and then reweigh it after withdrawing 1.0 mL of the amine from the reagent bottle. When this is done, the contents should weigh about 0.992 g.

Do the following steps with the help of the instructor or a trained assistant (see notes to the instructor) and the information in the Caution Box.

CAUTION

Butyllithium is pyrophoric and can burn in the presence of oxygen and moisture. The following operation must be conducted with an instructor or trained assistant. Anyone dispensing butyllithium must work in a hood and wear a lab coat, gloves, and goggles. Always use plastic syringes with needles equipped with *Luer Lock fittings when transferring butyllithium*. Make sure that the needle is securely attached to the syringe so that there is no danger that the needle will disconnect from the syringe. You must avoid any spills with butyllithium to avoid a fire. A dry chemical fire extinguisher should be immediately available in case of fires. The University of Minnesota has a guide for handling pyrophoric material. See pages 11–16 on the following website: http://www.dehs.umn.edu/PDFs/Pyrophoric_Chemicals_Guide.pdf

Obtain the butyllithium reagent bottle (1.6 M in hexanes) from the instructor.¹¹ Insert a *short syringe needle* from the inert gas source (gas manifold) into the reagent bottle equipped with a rubber septum. Position the needle so that it is above the liquid level in the butyllithium reagent bottle. This inert gas maintains a positive pressure and replaces the reagent withdrawn from the reagent bottle with the inert gas. Attach a *12-inch curved syringe needle with a Luer Lock fitting* to a 5-mL syringe also equipped with a Luer Lock adapter. Make sure that the needle is securely attached to the syringe so that there is no danger that the needle will disconnect from the syringe in the critical operations that follow. Puncture the septum on the butyllithium bottle with the syringe needle, but do not insert the needle into the solution yet. Keep the syringe in an inverted position (plunger down). Withdraw the inert gas from the bottle, and purge the needle/syringe three times with the inert gas, expelling the gas into the hood between each removal of inert gas. After the third purge, push the needle down into the liquid and withdraw 3.7 mL of butyllithium (0.0059 mol, 1.6 M in hexanes) from the reagent bottle while keeping the syringe inverted. Withdraw the needle slightly so that the tip is located above the liquid in the bottle and then draw in some inert gas to create an inert gas “blanket” above the butyllithium in the inverted syringe. This “blanket” will protect your butyllithium from contact with air (remember butyllithium is pyrophoric) when the syringe is withdrawn from the bottle. Now, remove the needle from the bottle and rapidly insert the needle into the rubber septum on the flask. Inject the butyllithium quickly through the rubber septum into your 250-mL flask. You will observe bubbles of butane gas that form when butyllithium reacts with dicyclohexylamine (see eq 1). The butane gas will escape through the exit syringe needle shown in Figure 17.11B. Remove the syringe and needle used to dispense the butyllithium from the rubber septum.¹² Stir the reaction mixture at room temperature for about 10 minutes or until evolution of butane ceases (you may need to close the stopcock temporarily on the gas manifold in order to see if bubbles of butane are being formed).

¹¹ Purchased **fresh each term** from Acros, n-butyllithium in hexanes, 1.6 M, #101433-282. Store this reagent in a refrigerator. Allow the butyllithium to warm to room temperature before using. Do not purchase any butyllithium with concentrations above 1.6 M, as they tend to be more dangerous to use. Never remove the rubber seal on the top of the bottle!

¹² When you are finished dispensing butyllithium, the syringe and needle will need to be cleaned. To clean the syringe and needle, draw up some hexane into the syringe and dispense the contents into 2-propanol (isopropyl alcohol). Repeat several times. Hexane dissolves the residual butyllithium and 2-propanol reacts with the butyllithium to neutralize the material.

Preparation of the Lithium Enolate (eq 2)

Using a 1-mL syringe equipped with a needle, draw up 0.70 mL of *tert*-butyl propionate ($C_7H_{14}O_2$, MW 130 g/mole, density 0.865 g/mL) into the syringe (0.605 g, 0.00465 mol).^{13,14} Insert the needle from the syringe containing the *tert*-butyl propionate into the rubber septum and add it in 0.1-mL portions over about 2 minutes. At this point, the lithium enolate of *tert*-butyl propionate forms (see eq 2).

Addition of the Aryl Bromide and Catalyst (eq 3,4,5,6)

Weigh out 0.989 g of 2-bromo-6-methoxynaphthalene (0.0042 mol).¹⁵ This material is not air sensitive, so you can leave it out in the open while you obtain the catalyst. Working quickly, weigh out the air-sensitive catalyst on a piece of *folded* weighing paper on an analytical (four place) balance to obtain about 0.0045 to 0.0055 g of dibromobis(tri-*tert*-butylphosphine)dipalladium¹⁶ $\{[Pd(t-Bu)_3PBr]_2\}$. Your instructor may do the weighing operation for you. This material is air sensitive so work as quickly as possible. Try to keep the catalyst weight in that range, but the most important thing is to work quickly.¹⁷

The following operations must be conducted as rapidly as possible since the flask will be open briefly to the atmosphere! Close the stopcock on the gas manifold, which maintains an inert atmosphere in the 250-mL flask. Now, working quickly, stop the stirrer, remove the inert gas needle from the manifold and the vent needle from the rubber septum. Remove the rubber septum from the flask and insert a powder funnel. Pour the 2-bromo-6-methoxynaphthalene as quickly as possible into the 250-mL flask. If some of the bromo compound adheres to the funnel, tap it a few times to dislodge as much as possible, and remove the funnel. Ignore any remaining solid in the funnel. Add the catalyst from the weighing paper directly into the opened flask. Immediately, reattach the rubber septum to the flask and insert the needle from the manifold and the exit needle into the septum. Open the stopcock on the gas manifold to reestablish an inert atmosphere in the flask. Start the stirrer again, and continue to flush the flask with inert gas for 15 to 30 minutes.¹⁸

Reaction Period

Remove the inert gas needle and vent needle from the rubber septum after 15 to 30 minutes, but *do not remove* the rubber septum. Store the round-bottom flask on a cork ring in a place in the laboratory recommended by your instructor. Allow the reaction to proceed for at least 24 hours or until the next laboratory period. During that time the solution becomes cloudy (formation of lithium bromide), and the mixture retains a yellow color. Be sure to maintain anhydrous conditions during this time period. The mechanism for the reaction of the lithium enolate with the 2-bromo-6-methoxynaphthalene is shown in eqs 4, 5, and 6.

¹³ Aldrich Chemical Company, *tert*-butyl propionate, #254525.

¹⁴ It may be helpful, but is not required, to preweigh the syringe and then reweigh the syringe after drawing up 0.70 mL into the syringe. This process gives an actual amount taken, of about 0.611 g.

¹⁵ Purchased from Alfa-Aesar, 2-bromo-6-methoxynaphthalene, #AAA19450-18.

¹⁶ Purchased from Alfa Aesar, dibromobis(tri-*tert*-butylphosphine)dipalladium #AA44446-77, purchased in 0.1-g quantities. This catalyst may also be named as palladium(I) tri-*tert*-butylphosphine bromide, dimer by other suppliers, $\{[Pd(t-Bu)_3PBr]_2\}$.

¹⁷ After the class has obtained the catalyst, the instructor should provide some inert gas to the remaining catalyst in the small vial. Cap the vial, cover the top with Parafilm, and place the vial back into a desiccator or in a refrigerator.

¹⁸ The solution may turn orange followed by a color change to a bronze or reddish color. Color changes are often variable, and are not predictable in all cases.

Working Up the Reaction (2nd Day of Lab Experiment)

After the reaction period, remove the rubber septum. You no longer need to maintain anhydrous conditions. Add 5 mL of 3 M aqueous HCl. Continue stirring for about 10 minutes. The precipitate that forms is dicyclohexylamine hydrochloride.

Dicyclohexylamine Hydrochloride Removal

Set up a 50-mm Büchner funnel for vacuum filtration using a 125-mL filter flask. Filter the mixture through a piece of Whatman #2 filter paper (4.25 cm). The filtration process may be slow. You can push the solid down a bit with a bent spatula to remove as much of the liquid as possible. *Keep the filtrate*, and add another 5 mL of 3 M aqueous HCl to the filtrate in the filter flask. The second addition of HCl produces more solid, but not as much as what you obtained after the first addition of HCl. Set up another vacuum filtration apparatus using your partner's 125-mL filter flask and 50-mm Büchner funnel. Refilter the mixture through a fresh piece of filter paper. Discard the solids that you collected on the filter papers. *Remember you want to keep the filtrates. Combine the two batches of filtrates. The filtrate consists of the product dissolved in toluene and the aqueous layer.*

Pour the contents of the filter flasks into a 125-mL separatory funnel. The upper organic layer contains toluene and the *tert*-butyl ester of naproxen. Remove the lower aqueous layer from the separatory funnel and discard it. Pour the toluene layer remaining in the funnel from the top of the separatory funnel into an Erlenmeyer flask.

Add a small amount of anhydrous sodium sulfate to dry the organic layer (see Technique 12, Section 12.9B). Decant the dried toluene layer into a 100-mL round-bottom flask.

Remove the toluene on a rotary evaporator, if available, or ask your instructor for advice on how to remove the toluene (a rotary evaporator is strongly advised because it is difficult to evaporate toluene by other methods). You may need to heat the water bath surrounding the flask to remove the toluene, but be gentle during the first part of the removal to avoid bumping the toluene up into the trap! This may take some time, as toluene has a relatively high boiling point even at reduced pressure. Remove the toluene until the level of liquid in the flask does not seem to be changing. The volume of liquid remaining will be very small, about 1 to 3 mL of liquid at this point.

Hydrolysis of the *tert*-butyl Ester of Naproxen with Trifluoroacetic Acid to Yield Naproxen

Add 4 mL of trifluoroacetic acid (**careful, toxic**) to the crude *tert*-butyl ester of naproxen contained in the round-bottom flask. Add a stir bar and heat the contents with stirring on a hot plate in an open flask at 130°C for ½ hour (check the internal temperature of the mixture with a thermometer). After ½ hour, cool the flask to room temperature, and add 20 mL water and 20 mL methylene chloride (dichloromethane).

Transfer the combined organic and aqueous phases into a separatory funnel, using small portions of additional methylene chloride to aid the transfer and to rinse the round-bottom flask. Shake the funnel to extract the naproxen into the organic layer. This time, the organic layer will be the lower layer (remember, methylene chloride is more dense than water). Remove the lower methylene chloride layer and save it. Discard the upper aqueous layer. Return the organic layer to the separatory funnel, extract it with another 20-mL portion of water, and again save the lower methylene chloride layer. Discard the aqueous layer. We now have the free carboxylic acid, naproxen, in the organic layer along with some side products.

Pour the methylene chloride layer back into the separatory funnel. At this point, the organic layer contains naproxen and miscellaneous organic compounds that need to be separated from the naproxen. Keeping in mind that naproxen is a carboxylic acid, we will now convert it to the *water soluble sodium salt of naproxen* to remove it from the remaining organic compounds that are not acidic.¹⁹ To accomplish the extraction of the naproxen as the sodium salt, add 10 mL of 5% aqueous NaOH to the separatory funnel. Shake the separatory funnel. The layers will separate slowly; wait until they separate completely. Drain the lower methylene chloride organic layer. *Pour the upper aqueous layer into a beaker from the top of the separatory funnel and save it, as it contains the water soluble sodium salt of naproxen.* Reintroduce the organic layer into the separatory funnel and reextract it with another 10-mL portion of 5% aqueous sodium hydroxide. Again, wait for the layers to separate and remove the lower methylene chloride layer. The organic layer may be discarded unless your instructor wants you to determine the GC/MS of the side products and starting material.²⁰ Pour the upper aqueous layer from the separatory funnel into the same beaker with the first extract. *Be sure to save the two combined aqueous layers.*

Obtain 10 mL of 3M aqueous HCl. Start adding it dropwise to the basic aqueous solution that contains the sodium salt of naproxen. Initially, the solution is heavily basic, and the naproxen is present as the sodium salt. As HCl is added, there will be localized formation of the white solid, naproxen. In that region, we have “neutralized” the solution. Using a stirring rod, mix the material so that the white solid redissolves. Overall, the solution is still basic. In effect we are using the appearance of the white solid as an indication of the end point. Think titration. Keep adding the HCl, dropwise, with stirring until the entire solution turns white (the end point!). At that point, check the pH to see if the solution is now acidic, about pH = 2. You may not use all of the HCl solution. Cool the mixture in an ice bath for 15 minutes.

Collect the naproxen (now in the acid form) on a Büchner funnel, on Whatman #2 filter paper, under vacuum. Use 5 mL of ice water to help transfer all of the solid from the beaker to the Büchner funnel. Carefully lift out the filter paper with the naproxen. Put the filter paper on a large watch glass. Dry the naproxen in a 115°C oven for about ½ hour. Scrape the naproxen off the filter paper and *weigh it*. Break it up and crush the solid into as fine a powder as possible. Put the naproxen back into the oven and heat it for another ½ hour.²¹ Again, weigh the solid to see if it has reached a constant weight. *It is very important to remove all of the water from the naproxen, so make sure that it has reached a constant weight.* Determine the melting point (typically 147–151°C) and calculate the percentage yield. At the option of the instructor, determine the infrared spectrum and the ¹H NMR spectrum.²² Compare your NMR spectrum to the one shown in Figure 1. This spectrum has expansions drawn as insets on the spectrum. The peaks are assigned on the structure shown. Note that sometimes the small quartet for the methine proton at 3.87 ppm may

¹⁹ The organic layer contains small amounts of unreacted 2-bromo-6-methoxynaphthalene, 2-methoxynaphthalene, dicyclohexylamine, and a very small amount of unhydrolyzed *tert*-butyl ester of naproxen.

²⁰ The GC–MS analysis will show 2-bromo-6-methoxynaphthalene, 2-methoxynaphthalene, dicyclohexylamine, and a trace of *tert*-butyl naproxen.

²¹ Alternatively, you may leave the solid naproxen to dry until the next laboratory period. In this case, be sure to leave the solid in an open container. Squishy solid contains too much water. Water leads to an unsuccessful resolution step. Water must be removed.

²² NMR spectra determined at lower field strengths will not resolve the aromatic protons, but the singlet for methoxy, doublet for methyl, and quartet for the methine protons should be resolved.

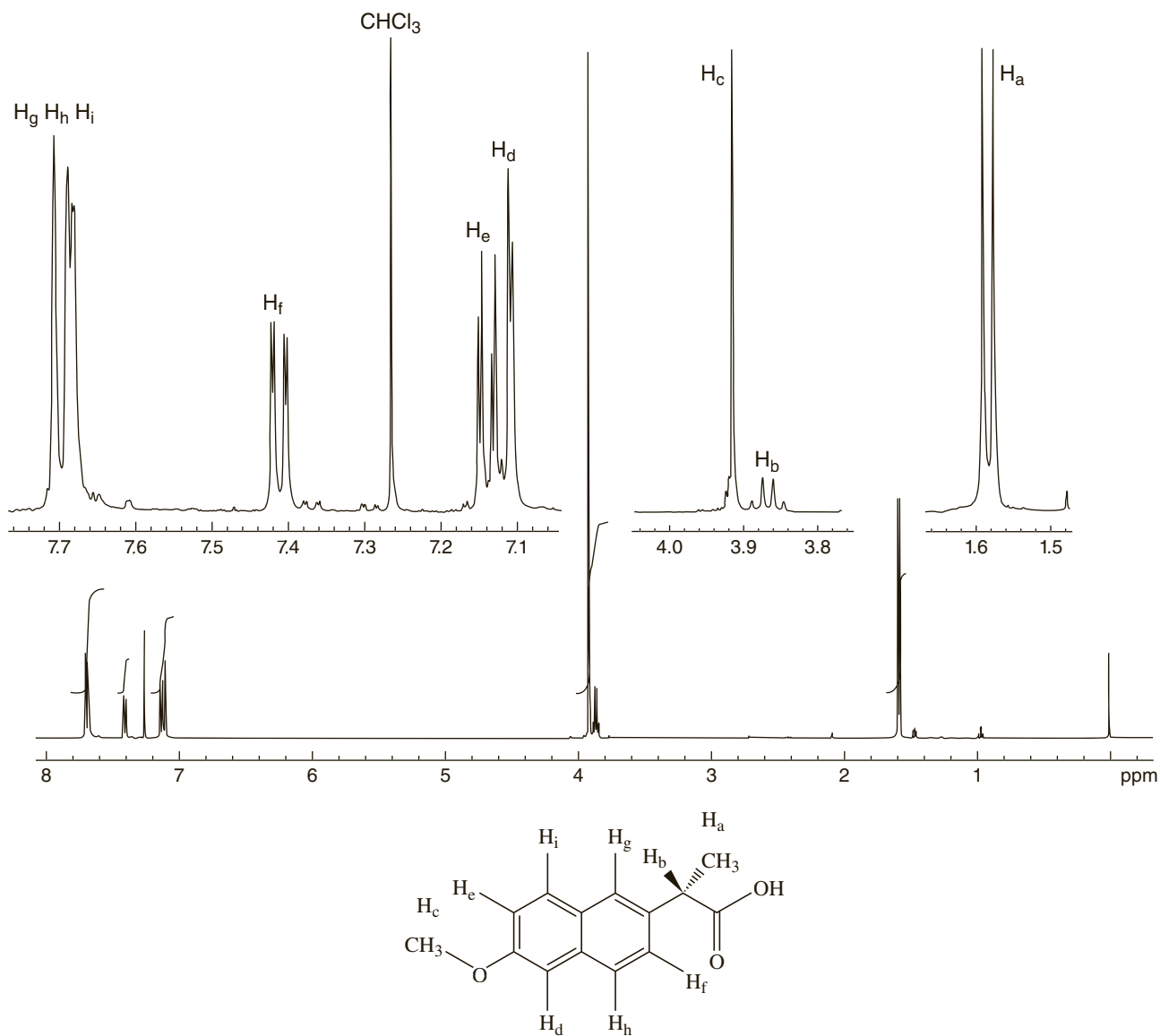


Figure 1


500 MHz ^1H NMR spectrum of racemic naproxen in CDCl_3 . The carboxylic acid proton does not appear in the spectrum. CHCl_3 appears at 7.26 ppm. Three regions in the spectrum have been expanded and shown as insets.

partially overlap with the large singlet for the methoxy group at 3.91 ppm. Also note that the carboxylic acid peak is not visible in the spectrum shown in Figure 1. During the same laboratory period, start the resolution of the naproxen to the (S)-enantiomer or wait until the next period. Save 0.2 g of racemic naproxen for the resolution step, and turn in the excess racemic naproxen to your instructor with your laboratory report.

Resolution of Racemic Naproxen: Isolation of the (S) Enantiomers (3rd Day of Lab Experiment)

This part of the experiment involves the resolution of the racemic naproxen to allow the isolation of the (S) enantiomer of naproxen. You will use the following procedure to

resolve the racemic naproxen. This procedure employs a chiral amine natural product, (*R*) (-)-cinchonidine,²³ to resolve racemic naproxen.²⁴


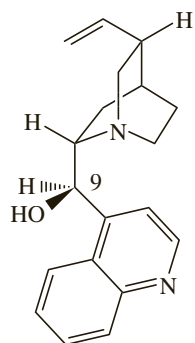


(-)-Cinchonidine, a naturally occurring alkaloid, from cinchona bark; has medicinal use as an antimalarial

Position 9 on the alkaloid has the (*R*) configuration.

MW = 294.4

In order to separate the enantiomers of naproxen, you need to form a pair of diastereomers. Diastereomers have different properties, including solubility characteristics.

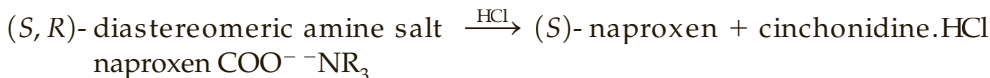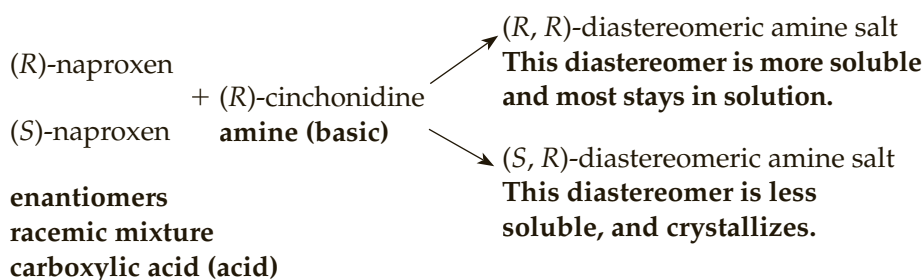

Procedure

Place 0.200 g of *dry racemic* naproxen and 0.261 g of (*R*)-(-) cinchonidine in a 25-mL Erlenmeyer flask. As shown above, two diastereomers will form. Add 6.0 mL of methanol and 2.5 mL of reagent grade acetone to the flask. Heat the contents of the flask slightly until all the solid dissolves. Now, stopper the flask, and allow the mixture to cool somewhat. Crystals of the less soluble (*S, R*) diastereomer should begin to form in a few minutes when the solution reaches room temperature. It is helpful to induce crystallization by scratching the inside of the flask with a stirring rod. One way to do this is to remove the stirring rod after scratching, blow on the end of the rod or wave the rod in the air to evaporate some of the solvent, and then reinsert the stirring rod back into the liquid. When you see a few crystals, *stopper the flask*, and allow the mixture to stand until the next laboratory period. If you do not see crystals, try seeding the solution with a few crystals from another student. Allow the mixture to stand until the next laboratory period. Try not to move or shake the flask, and do not cool the flask in an ice bath.

The next laboratory period (**4th day of lab experiment**), collect the (*S, R*) diastereomer crystals (some of the (*R, R*) diastereomer is also present) by vacuum filtration using a Büchner funnel or Hirsch funnel. Rinse the solid in the funnel with about 1 mL of ice-cold methanol. Dry the solid on the Büchner funnel or Hirsch funnel, under vacuum, for several minutes. Pour the remaining filtrate into a bottle that

²³ Purchased from Acros, cinchonidine, #200019-490.

²⁴ Harrison, I. T.; Lewis, B.; Nelson, P.; Rooks, W.; Roszkowski, A.; Tomolonis, A.; Fried, J. H. Nonsteroidal antiinflammatory agents. I. 6-Substituted 2-naphthylacetic acids. *Journal of Medicinal Chemistry*, **1970**, *13*(2), 203–205.

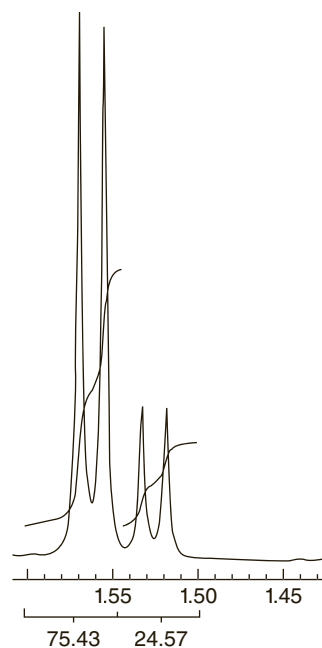


Figure 2

500 MHz ^1H NMR spectrum of a mixture of diastereomers, (*S*, *R*) and (*R*, *R*), formed from racemic naproxen and (*R*) cinchonidine. The expansion shows one doublet at 1.53 ppm for the methyl group in the (*R*, *R*) diastereomer and a larger doublet at 1.56 ppm for the methyl group in the (*S*, *R*) diastereomer. The integral values shown below the spectrum provide a value of 75.4% for the (*S*) enantiomer in the mixture.

your instructor will have made available to you. This filtrate contains mainly the (*R*, *R*) diastereomer.²⁵ Transfer the solid (*S*, *R*) diastereomer to a preweighed 25-mL Erlenmeyer flask and reweigh to determine the weight of the diastereomer formed from naproxen.²⁶ This material is enriched in the (*S*, *R*) diastereomer (amine salt), where (*S*) is the active enantiomer of naproxen and (*R*) is the enantiomer of cinchonidine. In principle, one could recrystallize this solid to further enrich the sample into even more of the (*S*, *R*) diastereomer, but this is a time-consuming process that will not be performed. *Be sure to save about 10 mg of this diastereomer for ^1H NMR analysis* (label the vial with your names and give it to the instructor). The 500 MHz ^1H NMR analysis of the diastereomers will yield the percentages of the (*S*) and (*R*) enantiomers present in the resolved sample of the (*S*, *R*) (major) and (*R*, *R*) (minor) diastereomers. The expansion shown in Figure 2 shows a pair of doublets for the methyl groups between 1.5 and 1.6 ppm for the two diastereomers. The taller doublet is assigned to the (*S*) enantiomer of naproxen (75.4%) while the smaller doublet is assigned to the (*R*) enantiomer of naproxen (24.6%). We will be able to compare these values to those obtained on the chiral HPLC column, described below.

Isolation of (*S*)-Naproxen

Add about 10 mL of methylene chloride to dissolve the remaining (*S*, *R*) diastereomer in the 25-mL Erlenmeyer flask. Add about 3 mL of water and add dropwise 1 mL of 6M aqueous HCl while swirling the contents of the flask. Using a stirring rod, remove

²⁵ The filtrate contains an enriched sample of the (*R*, *R*) diastereomer. Pour the filtrate into a bottle designated for this material. Your instructor may want to analyze the (*R*, *R*) diastereomer by NMR.

²⁶ You should expect about 0.120 to 0.240 g of the (*S*, *R*) diastereomer, starting from 0.200 g of racemic naproxen.

a drop of liquid from the flask and determine the pH to make sure that the solution is acidic. If the solution is still basic, add a few more drops of hydrochloric acid until the solution is acidic. Transfer the contents of the flask with a Pasteur pipette into a 60-mL separatory funnel. Rinse the flask with about 3 mL of water and transfer the water to the separatory funnel. Drain the lower, organic layer into a clean Erlenmeyer flask, and then reextract the aqueous layer in the separatory funnel with another 5-mL portion of methylene chloride. Combine the two methylene chloride extracts into the same flask. Add some (about 0.3 g) of anhydrous granular sodium sulfate to dry the methylene chloride extracts. Transfer the solution away from the drying agent, by carefully using a Pasteur pipette into a preweighed 50-mL round-bottom flask (your instructor may need to provide the flask). Use the rotary evaporator to remove the solvent. You now have an enriched sample of (*S*)-naproxen in the acid form (expect 26 to 77 mg). The (*S*)-naproxen will be analyzed on a chiral HPLC column. Dissolve 3 mg of your sample in 10 mL of the solvent provided (80:19.5:0.5 mixture of hexanes/isopropanol/acetic acid) in a snap-top container. Some of the (*R*) enantiomer will still be present in this sample. Your instructor may want to use the remaining resolved naproxen for polarimetry (see below). Turn in any excess resolved naproxen.

High-Performance Liquid Chromatography (HPLC)

You may use the following protocol for the separation of the enantiomers of naproxen with a chiral column on an HPLC instrument, if available. An (*S,S*) Whelk-O 1 (Regis Technologies, 4.6 mm × 250 mm, particle size: 5 microns) is equilibrated with an 80:19.5:0.5 mixture of hexanes/isopropanol/acetic acid at a flow rate of 2.0 mL/min on a Varian Pro Star chromatography HPLC system. Your instructor or assistant will manually inject a 25 μL portion of your resolved naproxen dissolved in an 80:19.5:0.5 mixture of hexanes/isopropanol/acetic acid into the HPLC system. The instrument is set up with a 2.0 mL/min flow rate, the elution is run for 10 min, and the UV peak is detected at 254 nm. An excellent separation of the two enantiomers is achieved. The (*R*)-(–)-naproxen elutes at 4.5 min and the (*S*)-(+)-naproxen elutes at 7.2 min.²⁷ In addition, the peak height at this concentration is about 0.15 absorbance units. Integration yields the percent composition of the enantiomers in the mixture. Your instructor will provide the HPLC chromatogram to you.²⁸ Use these percentages to determine the enantiometric excess (%ee) of the (*S*) enantiomer by the calculation shown in the first equation in Technique 23, Section 23.5 (replace moles with percentages in that equation). Attach the chromatogram to your report. Compare the percentages of the two enantiomers of naproxen obtained on the chiral HPLC column with the results for the percentages obtained on the mixture of resolved diastereomers obtained by NMR.

By the end of the lab period,

1. Be sure to submit a sample of the (*S*) enantiomer dissolved in 10 mL of the solvent (80:19.5:0.5 mixture of hexanes/isopropanol/acetic acid) in a snap-top container. This will be analyzed by chiral HPLC.
2. Also be sure to submit a 10 mg sample of the enriched (*S*, *R*) diastereomer for NMR analysis. Some of the (*R*, *R*) diastereomer will be present in the mixture.
3. Turn in your excess racemic naproxen in a labeled vial.
4. Turn in your excess resolved (*S*)-naproxen.

²⁷ You may also obtain retention times of anywhere from 3.7 to 4.1 minutes for the (*R*) enantiomer and from 5.4 to 6.1 minutes for the (*S*) enantiomer depending on the conditions of the chiral column.

²⁸ Typical chiral HPLC column results: 79 to 82% of the (*S*)-enantiomer.

Polarimetry (Optional)


The instructor may ask you to combine your remaining resolved naproxen with other students for determining the rotation of your (*S*)-naproxen by polarimetry. If so, your instructor will supply instructions. (*S*)-Naproxen has an observed specific rotation of $+66^\circ$. The solvent, chloroform, will be used as the solvent, unless you are told otherwise. Calculate the % optical purity (% enantiomeric excess) for your sample and compare the results with the chiral HPLC results.²⁹ Remember that the sample may only contain about 82% of the (*S*) enantiomers (Technique 23, Section 23.5) so you will not obtain a value of $+66^\circ$ from the polarimeter.

LABORATORY REPORT

1. Weight of the racemic naproxen, and percentage yield based on the 2-bromo-6-methoxynaphthalene.
2. Melting point of the racemic naproxen.
3. Interpretation of the infrared spectrum of racemic naproxen, if determined.
4. Comparison of the NMR spectrum of racemic naproxen with Figure 1, if determined.
5. Weight of the (*S, R*) diastereomer obtained in the resolution of racemic naproxen with cinchonidine.
6. Chiral HPLC results on naproxen, including percentages of each enantiomer.
7. Calculation of the enantiomeric excess of the (*S*)-enantiomer.
8. Report the percentages of each of the enantiomers in the NMR spectrum of the pair of diastereomers, (*S, R*) and (*R, R*), formed from (*R*)-cinchonidine. The (*S, R*) diastereomer is present in excess. Compare your values to those shown in Figure 2.

QUESTIONS

1. Would you expect the ^1H NMR spectrum of (*S*)-naproxen to be different from or identical with the NMR spectrum of racemic naproxen shown in Figure 1?
2. Give a balanced equation for the reaction of naproxen with sodium hydroxide to prepare sodium naproxen.
3. How could you prepare a sample of pure (*S*)-naproxen from a sample of racemic naproxen using chromatography?
4. How would the pattern differ for the methyl group patterns in the ^1H NMR spectrum of the (*R, R*)-diastereomer compared to that shown in Figure 2?
5. Draw the structure of the (*R*)-enantiomer of naproxen.
6. Draw the structure of dicyclohexylamine hydrochloride.
7. Draw the structure of *tert*-butyl ester of naproxen.
8. Show how you could synthesize racemic ibuprofen starting with 1-bromo-4-isobutylbenzene.

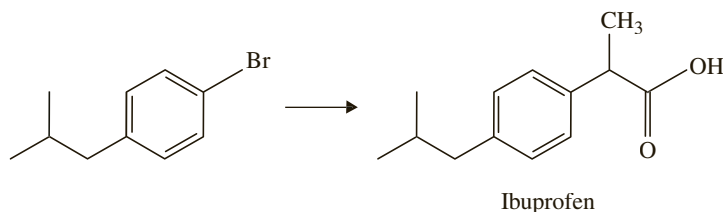

²⁹ Typical polarimetry results: 82% of the (*S*)-enantiomer. If the initially formed diastereomer is recrystallized, one may expect about 93% of the (*S*)-enantiomer.

The Aldehyde Enigma

Aldehyde chemistry

Extraction

Crystallization

Spectroscopy

Devising a procedure

Critical thinking application

The reaction mixture in this experiment contains 4-chlorobenzaldehyde, methanol, and aqueous potassium hydroxide. A reaction occurs that produces two organic compounds, Compound 1 and Compound 2. Both are solids at room temperature. Your task is to isolate, purify, and identify both compounds. A specific procedure is given for preparing the compounds, but you will need to work out the procedures for most other parts of this experiment.

SPECIAL INSTRUCTIONS

If the work on this experiment is done in pairs, work closely together as a team, dividing up the work equitably. A logical division of labor is for one student to work on Compound 1 and the other to work on Compound 2. Whether you work in pairs or not, you will need to plan your work carefully before coming to the laboratory, to make efficient use of class time.

SUGGESTED WASTE DISPOSAL

Dispose of all filtrates into the container designated for halogenated organic wastes.

PROCEDURE

This procedure should produce enough of each compound to complete the experiment; however, in some cases it may be necessary to run the reaction a second time. Although this experiment can be done individually, it works out especially well for two students to work together.

CAUTION

Be sure there is no acetone present on any of the glassware. Acetone will interfere with the desired reaction.

Running the Reaction

Add 1.50 g of 4-chlorobenzaldehyde and 4.0 mL of methanol to a 25-mL round-bottom flask. With gentle swirling, add 4.0 mL of an aqueous potassium hydroxide solution¹ with a Pasteur pipette. **Avoid getting potassium hydroxide solution on the ground-glass joint!** Add a stir bar to the flask and attach a water-cooled condenser. Using a hot water bath, heat the reaction mixture at about 65°C with stirring for 1 hour. Cool the mixture to room temperature and add 10 mL of water to the flask. Pour the mixture into a beaker and use another 10 mL of water to aid the transfer into the beaker.

Using a separatory funnel, extract the reaction mixture with 10 mL of methylene chloride. Drain the organic (bottom) layer into another container. Extract the aqueous layer with another 10-mL portion of methylene chloride. Combine the organic layers. The organic layer contains Compound 1, and the aqueous layer contains Compound 2.

Organic Layer

Wash the organic layer two times with 10-mL portions of 5% aqueous sodium bicarbonate solution. Then wash the organic layer with an equal volume of water. If an emulsion forms, use a little saturated sodium chloride solution to break it. Dry the organic layer over anhydrous sodium sulfate for 10–15 minutes. After the dried solution is removed from the drying agent, the organic layer should contain only Compound 1 and methylene chloride. Isolate Compound 1 by removing the methylene chloride.

Purify Compound 1 by crystallization. See “Testing Solvents for Crystallization,” Technique 11, Section 11.6, for instructions on how to determine an appropriate solvent. You should try 95% ethanol and xylene. After determining the best solvent, crystallize the compound using a hot water bath at about 70°C for heating to avoid melting the solid. Identify Compound 1 using some or all of the techniques given next in the section “Identification of Compounds.”

Aqueous Layer

To precipitate Compound 2, add 10 mL of cold water and acidify with 6M HCl. As acid is added, stir the mixture. Do not over-acidify the solution; pH 3 or 4 is fine. If no precipitate is formed on acidification, add saturated NaCl to aid the process. This is called **salting out**.

Isolate Compound 2 and dry it in an oven at about 110°C. Purify it by crystallization (see Technique 11, Section 11.5 for instructions on how to determine an appropriate solvent). You should try methanol and 95% ethanol. After determining the best solvent, purify the compound by crystallization and identify the purified solid using some or all of the techniques given in the next section, “Identification of Compounds.”

¹ Dissolve 61.7 g of potassium hydroxide in 100 mL of water.

IDENTIFICATION OF COMPOUNDS

Identify Compound 1 and Compound 2 using any of the following techniques:

1. *Melting point*: Consult a handbook for literature values.
2. *Infrared spectroscopy*: KBr pellet is preferred.
3. *Proton or carbon NMR*: Compound 1 dissolves easily in CDCl_3 ; use deuterated DMSO or Unisol to dissolve Compound 2.²
4. Some of the “wet” chemical tests listed in Experiment 57 may be helpful: solubility tests, Beilstein test for halide, and others you may think appropriate.
5. Physical properties such as color and shape of crystals may also be helpful.

REPORT

Write out a complete procedure by which you synthesized and isolated Compounds 1 and 2. Describe the results of your experiments to determine a good crystallization solvent for both compounds. Draw the structures of Compounds 1 and 2. Give all melting-point data and results of other tests used to identify the two compounds. Identify significant peaks in the infrared spectrum and proton/carbon NMR spectra. Show clearly how all these results confirm the identity of the two compounds. Write a balanced equation for the synthesis of Compounds 1 and 2. What type of reaction is this? Propose a mechanism for the reaction. Determine the percentage yield of each of the compounds.

² Unisol is a mixture of chloroform-d and DMSO-d₆ available from Norell, Inc., 120 Marlin Lane, Mays Landing, NJ 08330.

Synthesis of Substituted Chalcones: A Guided-Inquiry Experience

Crystallization


Aldol condensation

Use of the chemical literature

Project-based experiment

In Experiment 39, you were introduced to the **aldol condensation reaction**, which you used to prepare a variety of **benzalacetophenones** or **chalcones**. In this experiment, you will again prepare chalcones, but you will do so in a guided-inquiry experiment that simulates some of the methodology that you are likely to use in research.

You will select from a variety of substituted benzaldehydes (1) and substituted acetophenones (2) to prepare benzalacetophenones (chalcones) (3) that bear a combination of substituents in each of the aromatic rings (see figure).

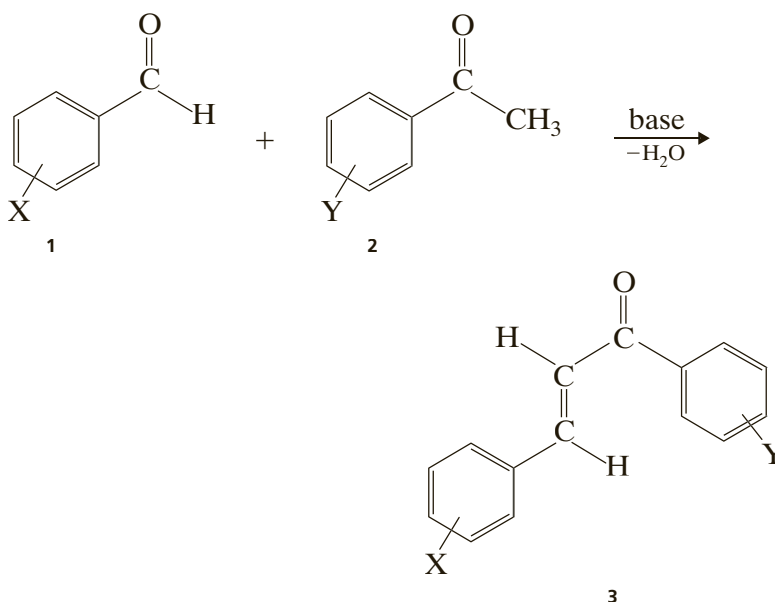

Once you have selected your starting materials, you will determine the complete structure of the condensation product that you expect to be formed in your reaction. You will also determine the molecular formula. With this information, you will be able to conduct an online literature search of *Chemical Abstracts* using **STN Easy** or **SciFinder Scholar**. From the literature search, you will be able to obtain the complete name of your target chalcone, its CAS Registry Number, and literature citations from the primary chemical literature. These literature citations should be

able to afford you characterization information about your target chalcone, including melting points, infrared spectra, and NMR spectra.

After you have conducted the literature search, the final step will be to prepare your chalcone and compare its properties with those that you were able to find in the literature.

The purpose of this experiment is to introduce you to many of the activities that you are likely to encounter in research. These include an examination of the target molecule, selection of appropriate starting materials, thorough search of the primary chemical literature, laboratory synthesis of the desired compound, and characterization (including a comparison of the physical properties of the product with published values found in journal articles or other tables of data).

REQUIRED READING

<i>Review:</i>	Technique 8	Filtration, Section 8.3
	Technique 11	Crystallization: Purification of Solids, Section 11.3
	Experiment 3	Crystallization
<i>New:</i>	Technique 29	Guide to the Chemical Literature

SPECIAL INSTRUCTIONS

Before beginning this experiment, you should select a substituted benzaldehyde and a substituted acetophenone. Your instructor will determine the method of assigning these reactants. You should also sign up for an STN Easy or SciFinder Scholar computer session. Your instructor will provide you with instructions on how to conduct an online computer search. Before coming to the computer session, you should work out the structure of your target compound and determine its molecular formula.

Note that sodium hydroxide solutions are caustic. Be careful when handling the substituted benzaldehydes and acetophenones. Wear personal protective equipment and work in a well-ventilated area.

SUGGESTED WASTE DISPOSAL


All filtrates should be poured into a waste container designated for nonhalogenated organic waste. Note that your instructor may establish a different method of collecting wastes for this experiment.

NOTES TO THE INSTRUCTOR


It is best to introduce this project two to three weeks before the date of the actual chalcone synthesis to allow time for searching the literature. You will have to develop a method of assigning a target compound to each student. You will also have to schedule computer time for the online searching of *Chemical Abstracts*. We recommend that you prepare a handout that describes how to search *Chemical Abstracts* using STN Easy or SciFinder Scholar. The handout should guide the students through the process of finding the registry number for the target compound and for finding pertinent references, with particular attention to

references describing the preparation of the compound. Finally, you will have to determine whether to require a formal laboratory report and what the expected format should be.

You may choose to create a multistep synthesis by linking this experiment to the Friedel-Crafts acylation reaction (Experiment 60) for the preparation of the substituted acetophenone. Experiment 60 contains suggestions for Friedel-Crafts acetophenones that work well when converted to chalcones. Following the synthesis of the chalcone in the current experiment, you can then carry out the cyclopropanation reaction (Experiment 67) and/or the epoxidation of the chalcone (Experiment 66). If the multistep scheme is to be followed, you should ask the class to scale up the chalcone preparation in order to have enough material to complete Experiment 66 and 67.

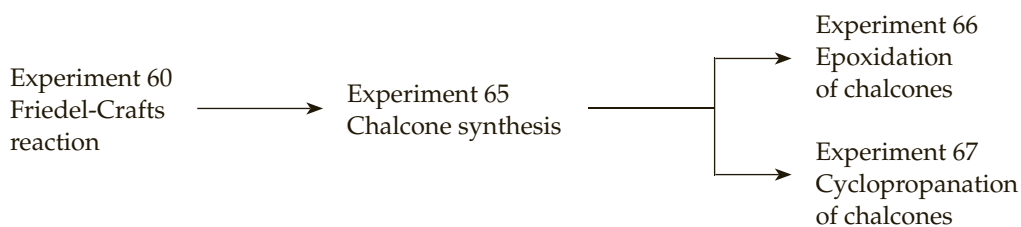

Another multistep synthesis, shown below, involves linking the synthesis of a chalcone in Experiment 65 with the epoxidation of the chalcone (Experiment 66) and/or the cyclopropanation of the chalcone (Experiment 67). If you plan for creating a multistep synthesis as described here, it may be a good idea to make a larger quantity of chalcone by scaling up the amounts of substituted acetophenone and substituted benzaldehyde used to prepare the chalcone.

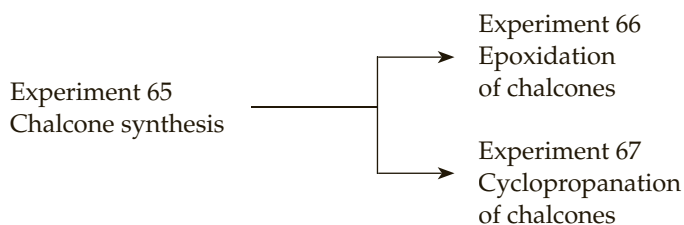

PROCEDURE

Before beginning the synthesis of your chalcone, determine its structure and molecular formula and perform the online search of *Chemical Abstracts*, following the instructions that your instructor provides.

Running the Reaction

Place 0.005 moles of your substituted benzaldehyde into a *tared* 50-mL Erlenmeyer flask, and reweigh the flask to determine the weight of material transferred.

Add 0.005 moles of the substituted acetophenone and 4.0 mL of 95% ethanol to the flask that contains the substituted benzaldehyde. Add a magnetic stirring bar to the flask. Swirl the flask to mix the reagents, and dissolve any solids present. It may be necessary to warm the mixture on a steam bath or hot plate to dissolve the solids. If this is necessary, the solution should be cooled to room temperature before proceeding to the next step.

Add 0.5 mL of sodium hydroxide solution to the benzaldehyde/acetophenone mixture.¹ Add a magnetic stir bar and stir the mixture. Before the mixture solidifies, you may observe some cloudiness. *Wait until the cloudiness has been replaced with an obvious precipitate settling out to the bottom of the flask before proceeding to the next paragraph.* Continue stirring until solid forms (approximately 3 to 5 minutes).² Scratching the inside of the flask with your microspatula or glass stirring rod may help to crystallize the chalcone.³

Isolation of the Crude Product

Add 10 mL of ice water to the flask *after a solid has formed as indicated in the previous paragraph.* Stir the solid in the mixture with a spatula to break up the solid mass. Transfer the mixture to a small beaker with 5 mL of ice water. Stir the precipitate to break it up, and then collect the solid, under vacuum, on a Hirsch or Büchner funnel. Wash the product with cold water. Allow the solid to air-dry for about 30 minutes. Weigh the solid.

Crystallization

Crystallize your entire sample from hot 95% ethanol. You will have to use the crystallization techniques introduced in Experiment 3 to crystallize the chalcone. Once the crystals have been allowed to dry thoroughly, weigh the solid, determine the percentage yield, and determine the melting point.

Spectroscopy

Determine the infrared spectrum of your product. Dissolve some of your chalcone in CDCl_3 (in some cases DMSO-d_6 may be required for sparingly soluble compounds) for ^1H NMR analysis. The chalcone spectrum will show a pair of doublets ($^3J = 16$ Hz appearing near 7.7 and 7.3 ppm) for the two vinyl protons in the starting chalcone. These vinyl protons in the chalcone appear in the same region as the aromatic protons on the benzene rings. However, the doublets for the protons in the benzene ring are more narrowly spaced ($^3J = 7$ Hz) than the doublets for the vinyl protons. Often you will see a singlet at 7.25 ppm for CHCl_3 present in the CDCl_3 solvent. In addition, a water peak may appear near 1.5 ppm. If deuterated DMSO had been used as solvent, you may see a pattern at about 2.6 ppm for non-deuterated DMSO. At the option of your instructor, determine the ^{13}C spectrum.

Laboratory Report


At the option of your instructor, you may be required to write a formal laboratory report. If this is the case, use the format that your instructor provides, or base your report on the style found in the *Journal of Organic Chemistry* (see Technique 29). If a literature search is required, use *SciFinder Scholar* to search for the melting point of your chalcone for comparison with the value you obtain. It should be noted here that when searching the chemical literature with *SciFinder Scholar*, you will find that *Chemical Abstracts* often does not use the name “chalcone” as the name of

¹ This reagent should be prepared in advance by the instructor in the ratio of 6.0 g of sodium hydroxide to 10 mL of water.

² In some cases, chalcone may not precipitate. If this is the case, stopper the flask and allow it to stand until the next laboratory period. It is sometimes helpful to add an additional portion of base. Usually chalcone will precipitate during that time.

³ In some cases, the aldol intermediate does not eliminate to form chalcone leading to an OH group in the infrared spectrum. In addition, chalcone may undergo a Michael addition of the enolate of the acetophenone on the chalcone. If either of these reactions occur, the ^1H NMR spectrum will show peaks in the 2.0–4.2 ppm range.

your compound. As an example, notice the name that is assigned to the following structure.

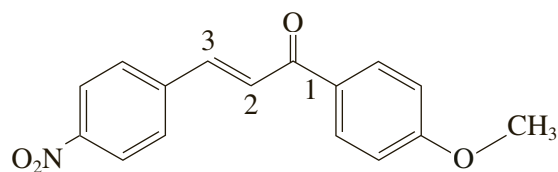

(E)-1-(4-methoxyphenyl)-3-(4-nitrophenyl)-2-propen-1-one

Submit the purified sample of your chalcone in a labeled vial to the instructor unless it is to be used in Experiment 66 and 67.


REFERENCES

Crouch, R. D.; Richardson, A.; Howard, J. L.; Harker, R. L.; and Barker, K. H. The Aldol Addition and Condensation: The Effect of Conditions on Reaction Pathway. *J. Chem. Educ.* **84** (2007): 475–476.

Vyvyan, J. R.; Pavia, D. L.; Lampman, G. M.; and Kriz, G. S. Preparing Students for Research: Synthesis of Substituted Chalcones as a Comprehensive Guided-Inquiry Experience. *J. Chem. Educ.* **79** (2002): 1119–1121.

QUESTIONS

1. Show how you begin with the indicated starting material and the Friedel-Crafts reaction to prepare the indicated chalcone products. You will require aldehydes and ketones in addition to the indicated starting material.

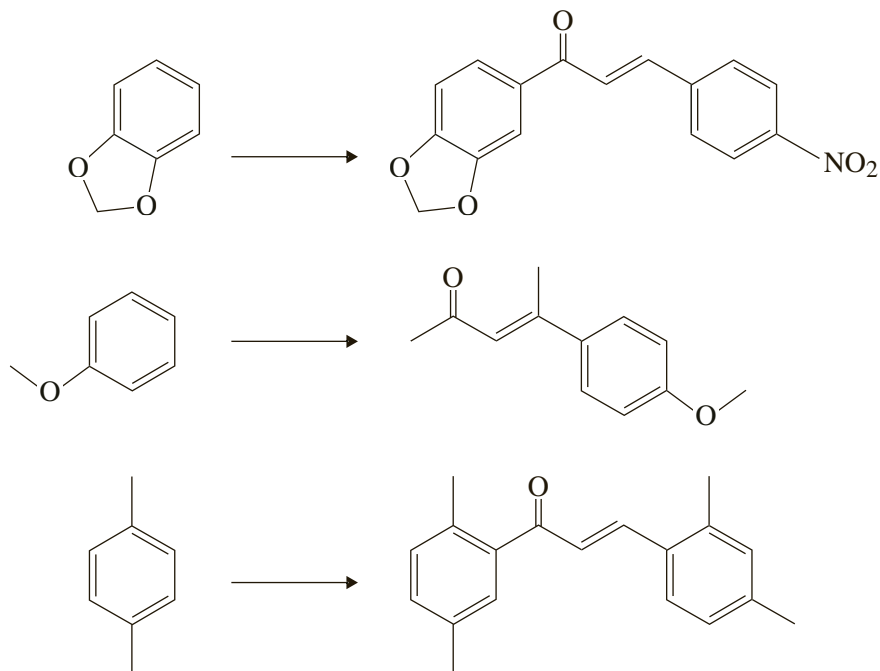


Green Epoxidation of Chalcones


Green chemistry

Reactions of chalcones

Epoxides are important intermediates widely used in multistep synthesis, and you have seen or will see them being used as intermediates in organic synthesis in your organic chemistry lecture courses. The common epoxidation reagent, *m*-chloroperoxybenzoic acid (*m*-CPBA), that you may have learned about in your lecture courses, does not work well on electron-poor conjugated ketones such as the chalcones employed in this experiment. Instead, we will use hydrogen peroxide in aqueous sodium hydroxide to prepare the epoxide. A “green” epoxidation of chalcones using these reagents has been reported in the literature, and this technique will be employed in this experiment.¹ The reaction is conducted in a not-so-green mixture of methanol, water, and dimethylsulfoxide (DMSO). DMSO is required to improve the solubility of the highly polar chalcones. The reaction mixture is stirred in an ice bath at 0°C for 1 hr to yield reasonable yields of epoxides. For example, *trans*-chalcone (1,3-diphenyl-2-propen-1-one) produces a 95% yield of the epoxide. Typical yields range from about 60 to 95% with other chalcones. To confirm that you have prepared the epoxide, you will analyze your product with ¹H NMR.

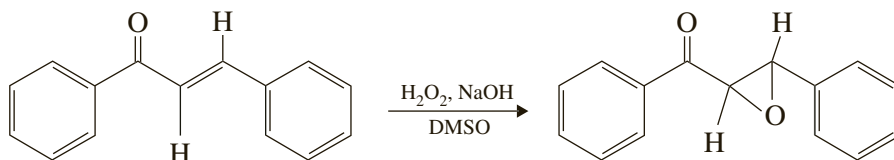

The mechanism follows the following pathway:

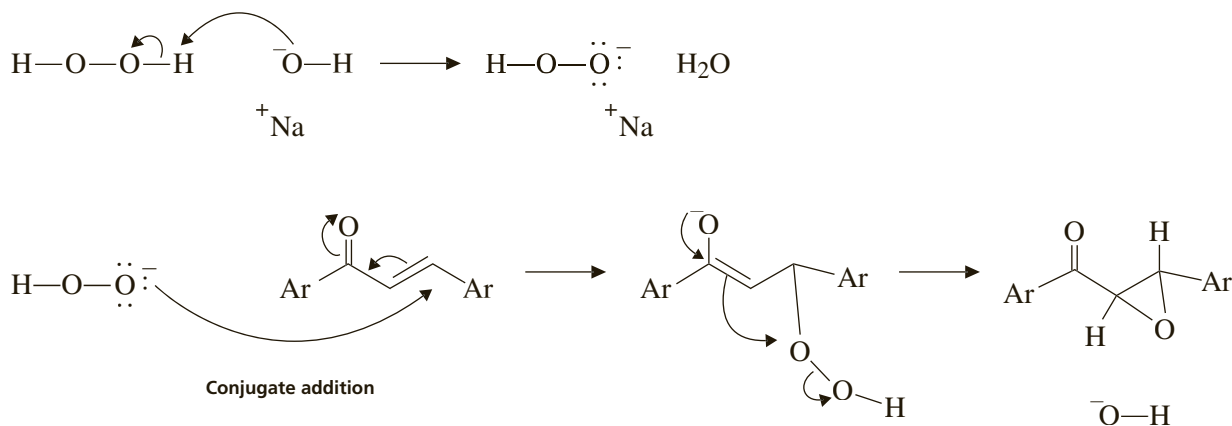

¹Fioroni, G.; Fringuelli, F.; Pizzo, F.; Vaccaro, L. Epoxidation of α, β -Unsaturated Ketones in Water. An Environmentally Benign Protocol. *Green Chemistry* **2003**, *5*, 425–428. Experiment developed by Butler, G., and Lampman, G.M., Western Washington University, Bellingham, WA.

REQUIRED READING

Review: Techniques 6, 7, 8, 11, 25, and 26

Read: Preparation of epoxides and reactions of epoxides in your lecture textbook

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous wastes in the container designated for aqueous waste. Place the organic waste in the nonhalogenated organic waste container.

NOTES TO THE INSTRUCTOR

Strong electron-releasing groups such as methoxy and methylenedioxy tend to retard the epoxide formation reaction on chalcones, leading to some residual chalcone remaining in the product. Alkyl groups are also electron-releasing, and they retard the formation of the epoxide. However, electron-withdrawing groups such as nitro and halogens enhance the reactivity of the chalcone. When halogen atoms are present along with methoxy, methylenedioxy, or alkyl groups, most of the chalcone is converted to the epoxide.

Students can determine the percent conversion of the chalcone to the epoxide by integrating one of the vinyl protons remaining in the aromatic region for the chalcone starting material and comparing that integral with the integral value for one of the protons on the epoxide ring. See the Spectroscopy section below for details.

PROCEDURE

Starting the Reaction

Add 0.50 mmole of your selected chalcone from Experiment 65, 3.5 mL of methanol, and a stir bar to a 50-mL round-bottom flask. Stir and gently heat the mixture for a few minutes to see if the chalcone will dissolve in methanol. If the chalcone does dissolve, proceed to the next paragraph. Most chalcones require some dimethylsulfoxide (DMSO) in addition to methanol to dissolve. Gradually add DMSO in 0.5-mL portions using a plastic Pasteur pipette until the chalcone dissolves with slight heating and stirring. It may take as much as 1 to 3 mL of DMSO to completely dissolve the chalcone. Now cool the round-bottom flask to room temperature. Some of the chalcone may precipitate as the temperature is lowered to room temperature, but the majority of the chalcone will remain in solution. You may proceed with the next step even if some solid remains.

Add 0.25 mL of 2 M aqueous sodium hydroxide using a plastic disposable pipette. Now add 65 μ L of 30% hydrogen peroxide using an automatic pipette. Support the flask in an ice-water bath with a clamp, but do not stopper the flask. Stir the mixture in an ice bath for 1 hr. Add more ice when necessary to keep the mixture between 0 and 2°C. Some chalcone will precipitate when the flask is cooled in the ice bath, but this should not be of concern because the chalcone will be converted to the epoxide even if some solid remains. Do not add any more DMSO.

Extraction with Diethyl Ether

Following the 1-hr reaction period, discontinue the stirring and add 5 mL of ice-cold water. A solid or possibly an oil should form. To extract the epoxide from the aqueous layer, add 10 mL of diethyl ether to the flask. Swirl the flask to help the epoxide

dissolve in diethyl ether. If necessary, add more diethyl ether to help dissolve the epoxide. The idea is to create two relatively clear layers, one aqueous and one organic. The amount of diethyl ether added is not critical.

Carefully transfer the mixture from the round-bottom flask to a separatory funnel. When pouring from the flask, use a funnel and a stir rod to direct the liquid into the funnel so that the liquid ends up in the funnel rather than on the surface of your hood! (It is difficult to pour from a flask with no lip!) Shake the funnel vigorously to extract the mixture, remove the lower aqueous layer, and pour the remaining ether layer from the top of the funnel into an Erlenmeyer flask. Now reintroduce the aqueous layer back into the separatory funnel, and re-extract it with another 10-mL portion of diethyl ether. After shaking, remove the lower aqueous layer, and again pour the ether extract from the top of the separatory funnel into the Erlenmeyer flask containing the first ether extract.

Drying and Removal of Diethyl Ether

Add anhydrous magnesium sulfate to the Erlenmeyer flask to dry the ether extracts. Cork the flask, and occasionally swirl the flask over a 5-min to 10-min period to dry the solution. Gravity-filter the solution through a piece of fluted filter paper into a preweighed 50- or 100-mL round-bottom flask (instructor-provided, if necessary). Remove the ether on the rotary evaporator, under vacuum. If a rotary evaporator is not available, your instructor will recommend an alternate method of removing solvent. A solid or an oil will form when the ether is removed. After the ether is removed, use a vacuum pump to remove the remaining solvent.

Isolation of the Epoxide

Reweigh the flask to determine the yield of the epoxide. Ideally, the isolated epoxide should be a solid, but often you will isolate an oily semi-solid (in the case of the oily semi-solid, proceed to the next paragraph). If a good-quality solid is obtained (check with your instructor for advice), add 8 mL of water to the solid to remove the DMSO that may have been extracted into ether. Bend the larger of the two spatulas you have in your drawer, and try to remove as much solid as possible from the sides and bottom of the round-bottom flask. Pour the solution containing the solid into a Hirsch or Büchner funnel attached to a filter flask, under vacuum, to collect the solid epoxide on filter paper. You may use additional cold water to aid the transfer process. Allow the solid to dry in an open container. When it is dry, weigh the solid and calculate the percentage yield. Also determine the melting point.

If the epoxide is an oily semi-solid, it will not be possible to collect the material on a Hirsch or Büchner funnel. Weigh the material and calculate the percentage yield. Dissolve the sample in CDCl_3 , and obtain the ^1H NMR spectrum as described in the next section.

Spectroscopy

Determine the infrared spectrum of your product. Dissolve some of your epoxide in CDCl_3 for ^1H NMR analysis. Compare the ^1H NMR spectrum of the starting chalcone with the spectrum of the epoxide. The starting chalcone spectrum will show a pair of doublets ($^3J = 16$ Hz appearing near 7.7 and 7.3 ppm) for the two vinyl protons in the starting chalcone. These vinyl protons in the chalcone appear in the same region as the aromatic protons on the benzene rings. However, the doublets for the protons in the benzene ring are more narrowly spaced ($^3J = 7$ Hz) than the doublets for the vinyl protons. The vinyl protons in the starting chalcone will be replaced with two peaks (actually a pair of doublets when expanded) near 4.0 to 4.4 ppm. The protons on the epoxide ring look like singlets in the NMR, but they are actually two finely spaced doublets ($^3J = 1.5$ to 2 Hz). Remember that you may see peaks in the spectrum for any remaining

DMSO at about 2.6 ppm. In addition, it is common to see a singlet for water appearing at about 1.5 ppm. At the option of your instructor, determine the ^{13}C spectrum.





Determine the percent conversion of the chalcone to the epoxide by integrating one of the vinyl protons that remains in the aromatic region for the chalcone starting material and comparing that integral with the integral value for one of the protons on the epoxide ring.

REFERENCES

- Dixon, C. E.; Pyne, S. G. Synthesis of Epoxidated Chalcone Derivatives. *J. Chem. Educ.* **1992**, *69*, 1032–1033.
- Fiorini, G.; Fringuelli, F.; Pizzo, F.; Vaccaro, L. Epoxidation of α , β -Unsaturated Ketones in Water: An Environmentally Benign Protocol. *Green Chem.* **2003**, *5*, 425–428.
- Fraile, J. M.; Garcia, J. I.; Mayoral, J. A.; Sebt, S.; Tahir, R. Modified Natural Phosphates: Easily Accessible Basic Catalyst for the Epoxidation of Electron-Deficient Alkenes. *Green Chem.* **2001**, *3*, 271–274.
- Maloney, G. P. Synthesis of 3-(2'-methoxy, 5'-bromophenyl)-2, 3-epoxyphenyl Propanone, a Novel Epoxidated Chalcone Derivative. *J. Chem. Educ.* **1990**, *67*, 617–618.


QUESTIONS

- Summarize the changes you expected to observe in the IR and ^1H NMR spectra of your epoxide product relative to the chalcone starting material.
- Draw the structures of the products expected in the following reactions.

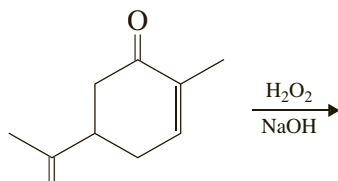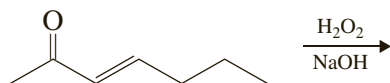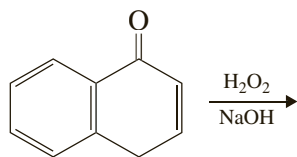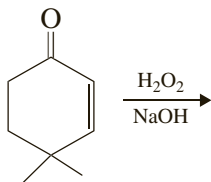

There are two $\text{C}=\text{C}$ double bonds, but only one reacts. Why?


- Draw the structure of the product expected in the following reaction.

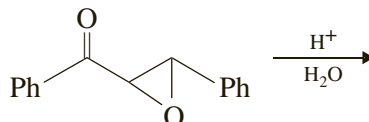


Cyclopropanation of Chalcones

Reaction of chalcones

The Corey–Chaykovsky reaction will be used to cyclopropanate your chalcone from Experiment 65. The reaction involves the reaction of trimethylsulfoxonium iodide and potassium *tert*-butoxide in anhydrous dimethylsulfoxide (DMSO).¹ The reaction is stirred at room temperature for 1 hr. For example, *trans*-chalcone (1,3-diphenyl-2-propen-1-one) produces an 88% yield of the cyclopropanated product. You will analyze your product by ¹H NMR and infrared spectroscopy.

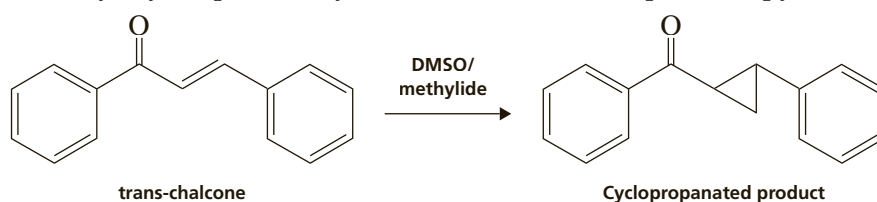

The mechanism follows the pathway below:

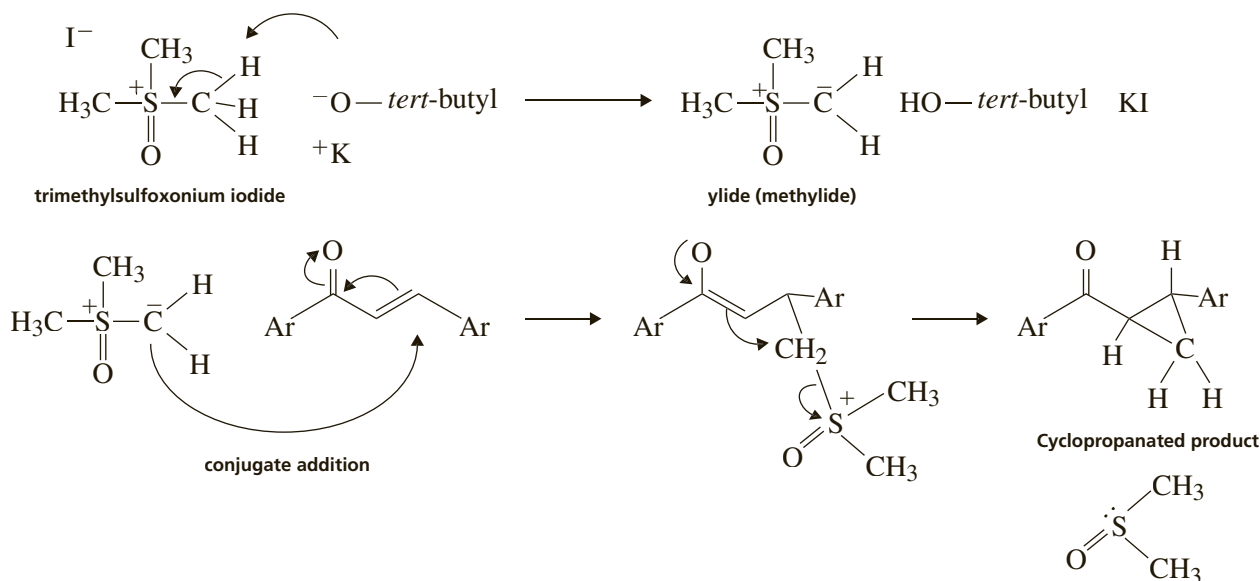

¹ Ciaccio, J. A.; Aman, C. E. Instant Methylyde Modification of the Corey–Chaykovsky Cyclopropanation Reaction. *Synthetic Communications* **2006**, *36*, 1333–1341. This experiment was developed by Truong, T. and Lampman, G. M., Western Washington University, Bellingham, WA.

REQUIRED READING

Review: Techniques 5, 6, 7, 8, 12, 20, 25, and 26

SUGGESTED WASTE DISPOSAL

Dispose of all aqueous wastes in the container designated for aqueous waste. Place the organic waste in the nonhalogenated organic waste container. Methylene chloride should be placed in the halogenated waste container.

NOTES TO THE INSTRUCTOR

Chalcones usually react completely in the cyclopropanation reaction, leaving little or no starting chalcone in the product.

PROCEDURE

Starting the Reaction

Dissolve the 0.50 mmole of the chalcone from Experiment 65 in 2.0 mL of *anhydrous* dimethylsulfoxide (DMSO)² in a 25-mL round-bottom flask. Allow the solid to dissolve.³ Add a stir bar. Add to the solution a dry mixture of Me₃S(O)I and KO-*tert*-butoxide (0.20 g, 0.6 mmol)⁴ in one batch. Now add a drying tube filled with CaCl₂ to the flask. Stir the solution for 1 hour at room temperature.

Extraction with Diethyl Ether

Transfer the mixture to a separatory funnel, and add 25 mL of saturated aqueous sodium chloride solution, using some of the sodium chloride solution to aid in the transfer of the reaction mixture to the funnel. Extract the mixture with a 15-mL portion of diethyl ether. Remove the lower aqueous layer, and pour the ether layer from the top of the separatory funnel into a beaker. Return the aqueous layer to the funnel, and re-extract it with another 15-mL portion of diethyl ether. Combine the two ether layers in the same beaker. Pour the ether extracts back into the separatory funnel, and re-extract the ether layer with two 25-mL portions of water, followed by extraction with 25 mL of saturated sodium chloride, each time draining the lower aqueous layer and saving the ether layer.

Drying and Removal of Diethyl Ether

Pour the diethyl ether layer from the top of the funnel into a dry Erlenmeyer flask, and dry the ether with anhydrous magnesium sulfate. Occasionally swirl the solution with the drying agent over a period of about 10 minutes. Gravity-filter the solution through a piece of fluted filter paper into a *preweighed* 50- or 100-mL

² Alfa Aesar, dimethyl sulfoxide, anhydrous, packaged under argon, Stock #43998, CAS #67-68-5

³ You may need to add more *anhydrous* DMSO to completely dissolve the chalcone.

⁴ The laboratory assistant should prepare the mixture by combining trimethylsulfoxonium iodide (Me₃S(O)I, 5.90 g; 26.8 mmol) with potassium *tert*-butoxide (KO-*tert*-Bu, 3.00 g; 26.7 mmol). Grind the mixture so that the two compounds are equally distributed and mixed with each other. One gram of this mixture provides 3.0 mmol of methylide/g or 0.6 mmole/0.2 g. Store the mixture in a desiccator.

round-bottom flask (instructor-provided, if necessary). Remove the ether on the rotary evaporator, under vacuum. If a rotary evaporator is not available, your instructor will recommend an alternative method of removing solvent. After the ether is removed, use a vacuum pump to remove any remaining ether from the sample. The product is likely to be an oil. Weigh the product, and determine the percentage yield.

Thin Layer Chromatography (Optional)

Check the purity of the product by TLC. Dissolve a small amount of the product in methylene chloride, and spot it on the plate. Also spot a dilute solution of the starting chalcone on the plate. Develop the plate in methylene chloride, and use the UV lamp to visualize the spots to see if there are any by-products or starting chalcone in your cyclopropanated product.

Spectroscopy

Determine the infrared spectrum of your product. Prepare an NMR sample for ^1H analysis in CDCl_3 . When the proton spectrum is returned to you, look for the disappearance of a pair of doublets ($^3J = 15$ Hz appearing near 7.7 and 7.3 ppm) for the vinyl protons in the starting chalcone (the normal expectation is for the chalcone to react completely). These doublets can be distinguished easily from the aromatic protons' doublets, which are more narrowly spaced ($^3J = 7$ Hz). These vinyl protons appear in the same region as the aromatic protons. The vinyl protons in the starting chalcone should be replaced by two cyclopropyl protons appearing at about 1.5 and 1.9 ppm for the diastereotopic protons in the CH_2 group. The two remaining cyclopropyl protons appear at about 2.6 and 2.88 ppm.⁵ At the option of your instructor, determine the ^{13}C spectrum.








REFERENCES

- Ciaccio, J. A.; Aman, C. E. Instant Methylide Modification of the Corey-Chaykovsky Cyclopropanation Reaction. *Synthetic Comm.* **2006**, *36*, 1333–1341.
- Corey, E. J.; Chaykovsky, M. Dimethyloxosulfonium Methylide and Dimethylsulfonium Methylide, Formation and Application to Organic Synthesis. *J. Am. Chem. Soc.* **1965**, *87*, 1353–1364.
- Lampman, G. M.; Koops, R. W.; Olden, C. C. Phosphorus and Sulfur Ylide Formation. *J. Chem. Educ.* **1985**, *62*, 267–268.
- Paxton, R. J.; Taylor, R. J. K. Improved Dimethylsulfoxonium Methylide Cyclopropanation Procedures, including a Tandem Oxidation Variant. *Synlett* **2007**, 633–637.
- Yanovskaya, L. A.; Dombrovsky, V. A.; Chizhov, O. S.; Zolotarev, B. M.; Subbotin, O. A.; Kucherov, V. F. Synthesis and Properties of *trans*-1-Aryl-2-Benzoylcyclopropanes and their Vinylogues. *Tetrahedron* **1972**, *28*, 1565–1573.

⁵ If instrumentation is available, run a gHSQC NMR experiment to confirm the assignment of the diastereotopic protons. This heteronuclear 2-D NMR experiment plots the carbon spectrum vs. the proton spectrum. The diastereotopic protons will correlate with only one ^{13}C peak at about 19 ppm. The other two cyclopropyl ring protons appear around 29 and 30 ppm in the ^{13}C spectrum.

QUESTIONS

- Summarize the changes you expect to observe in the IR and ^1H NMR spectra of your cyclopropane product relative to the chalcone starting material.
- Draw the structures of the products expected in the following reactions.

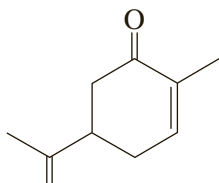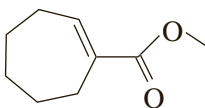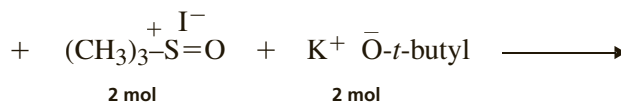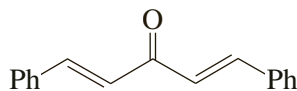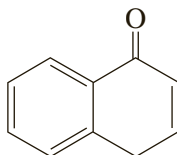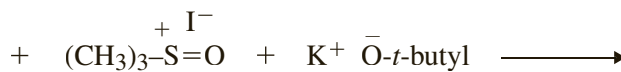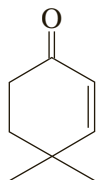

There are two $\text{C}=\text{C}$ double bonds, but only one reacts. Why?

Michael and Aldol Condensation Reactions

Aldol condensation


Michael reaction (conjugate addition)

Crystallization

Devising a procedure


Critical thinking application

In Experiment 39 (“The Aldol Condensation Reaction: Preparation of Benzalacetophenones”), substituted benzaldehydes are reacted with acetophenone in a crossed aldol condensation to prepare benzalacetophenones (chalcones). This is illustrated in the following reaction, where Ar and Ph are used as abbreviations for a substituted benzene ring and the phenyl group, respectively.

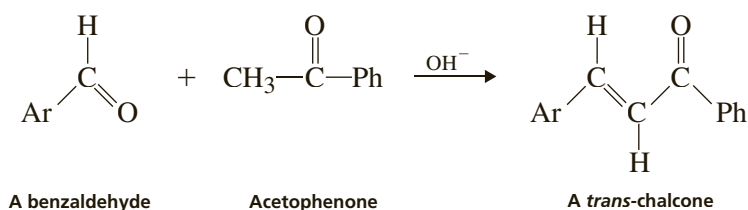

Experiment 41 involves the reaction between ethyl acetoacetate and *trans*-chalcone in the presence of base. Under the conditions of this experiment, a sequence of three reactions takes place: a Michael addition followed by an internal aldol reaction and a dehydration.

The purpose of this experiment is to combine the reactions introduced in Experiments 39 and 41 in the form of a project. Starting with one of four possible substituted benzaldehydes, you will synthesize a chalcone using the procedure given in Experiment 39. After performing a melting point to verify that this step has been completed successfully, you will perform a Michael/aldol reaction with the chalcone and ethyl acetoacetate using the procedure given in Experiment 41. The identity of this final product will be confirmed by melting point and possibly infrared and NMR spectroscopy.

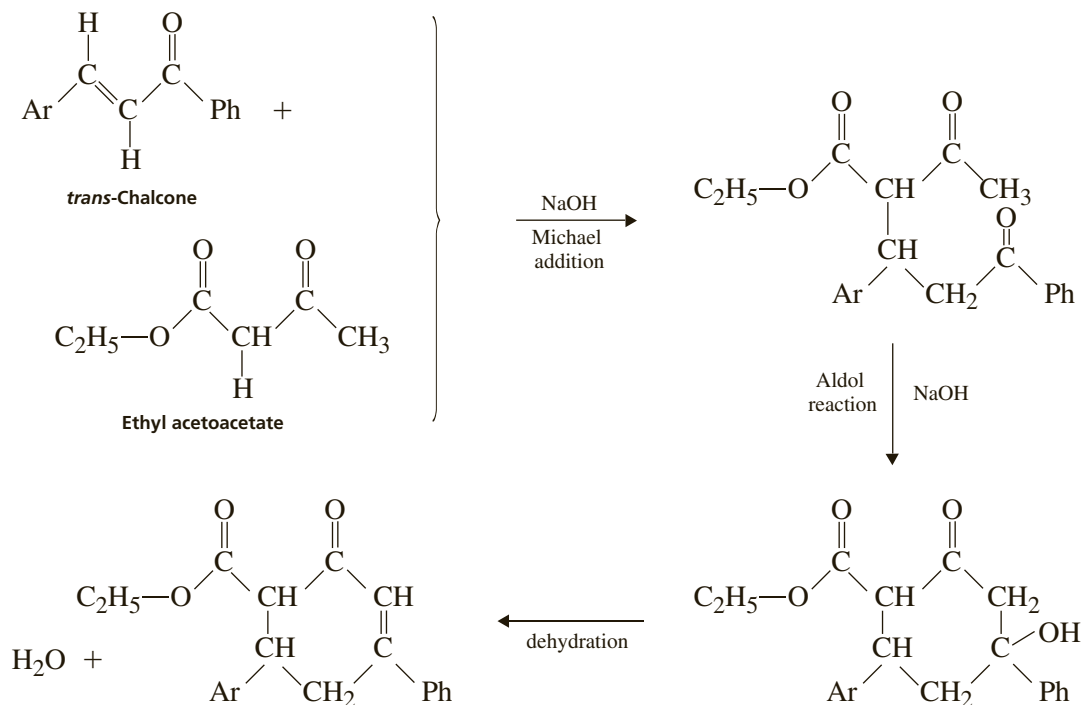

You will be assigned one of the aromatic aldehydes shown in the following list. For each aldehyde, the melting points of the corresponding chalcone and the Michael/aldol product are given:

Aldehyde	Chalcone (mp, °C)	Michael/Aldol Product (mp, °C)
4-Chlorobenzaldehyde	114–115	144–146
4-Methoxybenzaldehyde	73–74	107–109
4-Methylbenzaldehyde	95–97	150–152
Piperonaldehyde	119–121	151–153

REQUIRED READING

Review: Technique 11 Crystallization: Purification of Solids

SUGGESTED WASTE DISPOSAL

If your starting compound is 4-chlorobenzaldehyde, all filtrates should be poured into a waste container designated for halogenated organic wastes. If you use one of the other three aldehydes, dispose of all filtrates in the container designated for nonhalogenated organic wastes.

NOTES TO THE INSTRUCTOR

Some students may require individual help with this experiment. As a result, it may be difficult to use this experiment with a large class. It is a good idea to have

students prepare and present their procedure for approval before allowing them to begin the experimental work. The chalcones should be finely ground before being used in the second part of the experiment.

You may choose to have students react ethyl acetoacetate with one of the chalcones synthesized in Experiment 65. Because the product of the reaction may yield an unknown Michael/aldol product, a student will have the opportunity to conduct original research. A literature search may be incorporated with this exercise to see if the compound has been synthesized previously.

PROCEDURE

Your instructor will assign you one of the substituted benzaldehydes in the table above to use in this experiment. To prepare the chalcone, refer to the procedure in Experiment 39. To convert the chalcone to the Michael/aldol product, refer to the procedure given in Experiment 41. Using these procedures as a guide, devise the entire experimental procedure together with reagent quantities. The chalcone you prepare should be finely ground before using it in the second part of this experiment.

Initially, you should follow the procedures in Experiments 39 and 41 as closely as possible with appropriate adjustments in the scale. However, there is one part of the procedure in Experiment 41 that must be modified (see “Removal of Catalyst” in Experiment 41). The purpose of adding acetone in this step is to dissolve your product, leaving the solid catalyst behind. Depending on which substituted benzaldehyde you started with, different volumes of acetone may be required. Rather than following the instructions to add 1.5 mL of acetone, you should add a smaller portion and then stir with a spatula to see if most of the solid dissolves. If it does not, continue to add more acetone in small portions while stirring the mixture. When it is clear that most of the solid has dissolved and more solvent does not dissolve any more solid, then you can stop adding acetone. It is likely that you will need to add more than 1.5 mL of acetone, assuming the same scale as in Experiment 41.

If either procedure in Experiment 39 or 41 does not work, you may need to modify the procedure and run the experiment again. An unsuccessful procedure will most likely be indicated by either the melting point or spectral data. The problem you would most likely encounter in preparing the chalcone is difficulty in getting the product to solidify from the reaction mixture. The Michael/aldol reaction is more complicated, because there are two intermediate compounds that could be present in a significant amount in the final sample. If this occurs, both the melting point and the infrared spectrum may provide clues about what happened. It is possible you will need to increase the reaction time for this part of the experiment.

You must pay attention to scale so that you prepare enough of the chalcone for use in the next step and so that you finish up with a reasonable amount of the final product, about 0.1–0.2 g. It is possible, therefore, that the amounts of reagents given in Experiments 39 and 41 will need to be adjusted. If the scale needs to be changed in either experiment, be sure to adjust the amounts of all reagents proportionately and make any necessary changes in the glassware. In making your initial decision about scale, assume that the percentage yield of the chalcone after crystallization will be about 50%. Likewise, assume that the procedure in Experiment 41 will result in about a 50% yield.

To determine an accurate melting point of the chalcone or the final product, the sample must be pure and dry. In most cases, 95% ethanol can be used to crystallize these compounds. If this solvent does not work, you can use the procedure in

Technique 11, Section 11.6, to find an appropriate solvent. Other solvents to try include methanol or a mixture of ethanol and water. If you are unsuccessful in finding an appropriate solvent, consult your instructor.

It is particularly important that the chalcone be highly pure before going on to the next step. When you determine the amount of hot solvent to add when crystallizing the chalcone, it is best to add **more** than the minimum amount required to dissolve the solid. Otherwise, the amount of mother liquor may be so small that many of the impurities will not be removed during the vacuum filtration step. If the melting point after crystallization is not within 3–4°C of the melting point given in the table at the beginning of this experiment, you may need to crystallize the material a second time.

SPECTROSCOPY

Infrared Spectrum

You should obtain an infrared spectrum of the chalcone and the final product to verify the identity of each product in the reaction sequence. Obtain the infrared spectrum by the dry film method (Technique 25, Section 25.4) or as a KBr pellet (Technique 25, Section 25.5). For the Michael/aldol product, you should observe absorbances at about 1735 and 1660 cm^{-1} for the ester carbonyl and enone groups, respectively.

NMR Spectrum

Your instructor may also want you to determine the proton and carbon NMR spectra of each product. These may be run in CDCl_3 solvent. Some of the expected signals can be determined by referring to the NMR spectrum shown in Figure 2, Experiment 41. Although these data are for a slightly different compound, many of the signals will have similar splitting patterns and similar chemical shifts.

REPORT

The report should include balanced equations for the preparation of the chalcone and the Michael/aldol product. You should calculate both the theoretical and percentage yields for each step. Write your complete procedure as you actually performed it. Include the actual results of your melting-point determinations, and compare them to the expected results.

Include any infrared spectra obtained, and interpret the major absorption peaks. If you determined NMR spectra, you should include them, along with an interpretation of the peaks and splitting patterns.

REFERENCE

Garcia-Raso, A.; Garcia-Raso, J.; Campaner, B.; Maestres, R; Sinisterra, J. V. An Improved Procedure for the Michael Reaction of Chalcones. *Synthesis*. **1982**, 1037.

Esterification Reactions of Vanillin: The Use of NMR to Solve a Structure Proof Problem¹

Esterification


Crystallization

Use of a Craig tube

Nuclear magnetic resonance

Critical thinking application

The reaction of vanillin with acetic anhydride, in the presence of base, is an example of the esterification of a phenol. The product, which is a white solid, can be characterized easily by its infrared and NMR spectra.

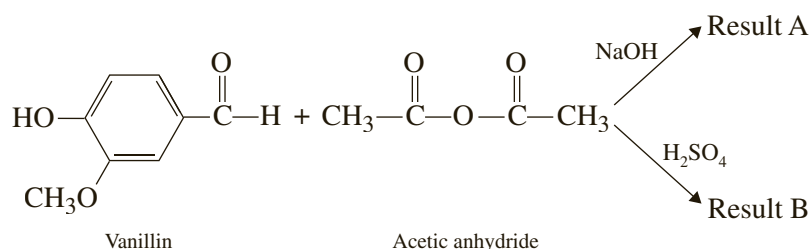

When vanillin is esterified with acetic anhydride under acidic conditions, however, the product that is isolated has a different melting point and different spectra. The object of this experiment is to identify the products formed in each of these reactions and to propose mechanisms that will explain why the reaction proceeds differently under basic and acidic conditions.

REQUIRED READING

Review: Techniques 8, 11, 25, and 26

You should also read the sections in your lecture textbook that deal with the formation of esters and nucleophilic addition reactions of aldehydes.

¹ This experiment is based on a paper presented at the 12th Biennial Conference on Chemical Education, Davis, California, August 2–7, 1992, by Professor Rosemary Fowler, Cottey College, Nevada, Missouri. The authors are very grateful to Professor Fowler for her generosity in sharing her ideas.

SPECIAL INSTRUCTIONS

Sulfuric acid is corrosive. Do not allow it to touch your skin.

SUGGESTED WASTE DISPOSAL

All filtrates and organic residues should be disposed of into the container designated for nonhalogenated organic wastes. Dispose of solutions used for NMR spectroscopy in the waste container designated for the disposal of halogenated materials.

PROCEDURE

Preparation of 4-Acetoxy-3-Methoxybenzaldehyde (Vanillyl Acetate)

Dissolve 0.30 g of vanillin in 5 mL of 10% sodium hydroxide in a 50-mL Erlenmeyer flask. Add 6 g of crushed ice and 0.8 mL of acetic anhydride. Stopper the flask with a clean, rubber stopper and shake it several times over a 20-minute period. A cloudy, milky white precipitate will form immediately on adding the acetic anhydride. Filter the precipitate, using a Hirsch funnel, and wash the solid with three 1-mL portions of ice-cold water.

In a Craig tube, recrystallize the solid from 95% ethyl alcohol. Heat the mixture in a hot water bath at about 60°C to avoid melting the solid. When the crystals are dry, weigh them and calculate the percentage yield. Obtain the melting point (literature value is 77–79°C). Determine the infrared spectrum of the product using the dry film method. Determine the proton NMR spectrum of the product in CDCl₃ solution. Using the spectral data, confirm that the structure of the product is consistent with the predicted result.

Esterification of Vanillin in the Presence of Acid

Place a magnetic spin vane in a 3-mL conical vial. Add 0.15 g of vanillin and 1.0 mL of acetic anhydride to the conical vial. Stir the mixture at room temperature until the solid dissolves. While continuing to stir the mixture, using a Pasteur pipette add 1 drop of 1.0M sulfuric acid to the reaction mixture. Cap the vial and stir at room temperature for 1 hour. During this period, the solution will turn purple or purple-orange.

At the end of the reaction period, transfer the reaction mixture to a centrifuge tube with a screw cap. Cool the tube in an ice-water bath for 3–4 minutes. Add 3.5 mL of ice-cold water to the mixture in the centrifuge tube. Cap the tube and shake it vigorously—almost as hard as you can shake! Continue to cool and shake the tube to induce crystallization. Crystallization has occurred when you can see small solid clumps separating from the cloudy liquid and settling to the bottom of the tube. If crystallization does not occur after 10–15 minutes, it may be necessary to seed the mixture with a small crystal of the product. Filter the product on a Hirsch funnel and wash the solid with three 1-mL portions of ice-cold water.

Using a Craig tube, recrystallize the crude product from hot 95% ethanol. Allow the crystals to dry. Weigh the dried crystals, calculate the percentage yield, and determine the melting point (literature value is 90–91°C). Determine the infrared spectrum of the product using the dry film method. Determine the proton NMR spectrum of the product in CDCl₃ solution.

REPORT

Compare the two sets of spectra obtained for the base- and acid-promoted reactions. Using the spectra, identify the structures of the compounds formed in each reaction. Record the melting points and compare them to the literature values. Write balanced equations for both reactions and calculate the percentage yields. Outline mechanistic pathways to account for the formation of both products isolated in this experiment.

The Techniques

Laboratory Safety

In any laboratory course, familiarity with the fundamentals of laboratory safety is critical. Any chemistry laboratory, particularly an organic chemistry laboratory, can be a dangerous place in which to work. Understanding potential hazards will serve you well in minimizing that danger. It is ultimately your responsibility, along with your laboratory instructor's, to make sure that all laboratory work is carried out in a safe manner.

1.1 Safety Guidelines It is vital that you take necessary precautions in the organic chemistry laboratory. Your laboratory instructor will advise you of specific rules for the laboratory in which you work. The following list of safety guidelines should be observed in all organic chemistry laboratories.

A. Eye Safety

Always Wear Approved Safety Glasses or Goggles. It is essential to wear eye protection whenever you are in the laboratory. Even if you are not actually carrying out an experiment, a person near you might have an accident that could endanger your eyes. Even dishwashing can be hazardous. We know of cases in which a person has been cleaning glassware—only to have an undetected piece of reactive material explode, throwing fragments into the person's eyes. To avoid such accidents, wear your safety glasses or goggles at all times.

Learn the Location of Eyewash Facilities. If there are eyewash fountains in your laboratory, determine which one is nearest to you before you start to work. If any chemical enters your eyes, go immediately to the eyewash fountain and flush your eyes and face with large amounts of water. If an eyewash fountain is not available, the laboratory will usually have at least one sink fitted with a piece of flexible hose. When the water is turned on, this hose can be aimed upward, and the water can be directed into the face, working much as an eyewash fountain does. To avoid damaging the eyes, the water flow rate should not be set too high, and the water temperature should be slightly warm.

B. Fires

Use Care with Open Flames in the Laboratory. Because an organic chemistry laboratory course deals with flammable organic solvents, the danger of fire is frequently present. Because of this danger, DO NOT SMOKE IN THE LABORATORY. Furthermore, use extreme caution when you light matches or use any open flame. Always check to see whether your neighbors on either side, across the bench, and behind you are using flammable solvents. If so, either wait or move to a safe location, such as a fume hood, to use your open flame. Many flammable organic substances are the source of dense vapors that can travel for some distance down a bench. These vapors present a fire danger, and you should be careful, as

the source of those vapors may be far away from you. Do not use the bench sinks to dispose of flammable solvents. If your bench has a trough running along it, pour only *water* (no flammable solvents!) into it. The troughs and sinks are designed to carry water—not flammable materials—from the condenser hoses and aspirators.

Learn the Location of Fire Extinguishers, Fire Showers, and Fire Blankets. For your own protection in case of a fire, you should immediately determine the location of the nearest fire extinguisher, fire shower, and fire blanket. You should learn how to operate these safety devices, particularly the fire extinguisher. Your instructor can demonstrate this.

If there is a fire, the best advice is to get away from it and let the instructor or laboratory assistant take care of it. **DON'T PANIC!** Time spent thinking before acting is never wasted. If it is a small fire in a container, it can usually be extinguished quickly by placing a wire-gauze screen with a ceramic fiber center or, possibly, a watch glass over the mouth of the container. It is good practice to have a wire screen or watch glass handy whenever you are using a flame. If this method does not extinguish the fire and if help from an experienced person is not readily available, then extinguish the fire yourself with a fire extinguisher.

Should your clothing catch on fire, **DO NOT RUN.** Walk *purposefully* toward the fire shower station or the nearest fire blanket. Running will fan the flames and intensify them.

C. Organic Solvents: Their Hazards

Avoid Contact with Organic Solvents. It is essential to remember that most organic solvents are flammable and will burn if they are exposed to an open flame or a match. Remember also that on repeated or excessive exposure, some organic solvents may be toxic, carcinogenic (cancer causing), or both. For example, many chlorocarbon solvents, when accumulated in the body, result in liver deterioration similar to cirrhosis caused by excessive use of ethanol. The body does not easily rid itself of chlorocarbons nor does it detoxify them; they build up over time and may cause future illness. Some chlorocarbons are also suspected of being carcinogens. **MINIMIZE YOUR EXPOSURE.** Long-term exposure to benzene may cause a form of leukemia. Do not sniff benzene and avoid spilling it on yourself. Many other solvents, such as chloroform and ether, are good anesthetics and will put you to sleep if you breathe too much of them. They subsequently cause nausea. Many of these solvents have a synergistic effect with ethanol, meaning that they enhance its effect. Pyridine causes temporary impotence. In other words, organic solvents are just as dangerous as corrosive chemicals, such as sulfuric acid, but manifest their hazardous nature in other, more subtle ways.

If you are pregnant, you may want to consider taking this course at a later time. Some exposure to organic fumes is inevitable, and any possible risk to an unborn baby should be avoided.

Minimize any direct exposure to solvents and treat them with respect. The laboratory room should be well ventilated. Normal cautious handling of solvents should not result in any health problems. If you are trying to evaporate a solution in an open container, you must do the evaporation in the hood. Excess solvents should be discarded in a container specifically intended for waste solvents, rather than down the drain at the laboratory bench.

A sensible precaution is to wear gloves when working with solvents. Gloves made from polyethylene are inexpensive and provide good protection. The disadvantage of polyethylene gloves is that they are slippery. Disposable surgical gloves provide a better grip on glassware and other equipment, but they do not offer as much protection as polyethylene gloves. Nitrile gloves offer better protection.

Do Not Breathe Solvent Vapors. In checking the odor of a substance, be careful not to inhale very much of the material. The technique for smelling flowers is not advisable here; you could inhale dangerous amounts of the compound. Rather, a technique for smelling minute amounts of a substance should be used. Pass a stopper or spatula moistened with the substance (if it is a liquid) under your nose. Or hold the substance away from you and waft the vapors toward you with your hand. But *never* hold your nose over the container and inhale deeply!

The hazards associated with organic solvents you are likely to encounter in the organic laboratory are discussed in detail in Section 1.3. If you use proper safety precautions, your exposure to harmful organic vapors will be minimized and should present no health risks.

Safe Transportation of Chemicals. When transporting chemicals from one location to another, particularly from one room to another, it is always best to use some form of **secondary containment**. This means that the bottle or flask is carried inside another, larger container. This outer container serves to contain the contents of the inner vessel in case a leak or breakage should occur. Scientific suppliers offer a variety of chemical-resistant carriers for this purpose.

D. Waste Disposal

Do Not Place Any Liquid or Solid Waste in Sinks; Use Appropriate Waste Containers. Many substances are toxic, flammable, and difficult to degrade; it is neither legal nor advisable to dispose of organic solvents or other liquid or solid reagents by pouring them down the sink.

The correct disposal method for wastes is to put them in appropriately labeled waste containers. These containers should be placed in the hoods in the laboratory. The waste containers will be disposed of safely by qualified persons using approved protocols.

Specific guidelines for disposing of waste will be determined by the people in charge of your particular laboratory and by local regulations. Two alternative systems for handling waste disposal are presented here. For each experiment that you are assigned, you will be instructed to dispose of all wastes according to the system that is in operation in your laboratory.

In one model of waste collection, a separate waste container for each experiment is placed in the laboratory. In some cases, more than one container, each labeled according to the type of waste that is anticipated, is set out. The containers will be labeled with a list that details each substance that is present in the container. In this model, it is common practice to use separate waste containers for aqueous solutions, organic halogenated solvents, and other organic nonhalogenated materials. At the end of the laboratory class period, the waste containers are transported to a central hazardous materials storage location. These wastes may be later consolidated and poured into large drums for shipping. Complete labeling, detailing each chemical contained in the waste, is required at each stage of this waste-handling process, even when the waste is consolidated into drums.

In a second model of waste collection, you will be instructed to dispose of all wastes in one of the following ways:

Nonhazardous solids. Nonhazardous solids such as paper and cork can be placed in an ordinary wastebasket.

Broken glassware. Broken glassware should be put into a container specifically designated for broken glassware.

Organic solids. Solid products that are not turned in or any other organic solids should be disposed of in the container designated for organic solids.

Inorganic solids. Solids such as alumina and silica gel should be put in a container specifically designated for them.

Nonhalogenated organic solvents. Organic solvents such as diethyl ether, hexane, and toluene, or any solvent that does not contain a halogen atom, should be disposed of in the container designated for nonhalogenated organic solvents.

Halogenated solvents. Methylene chloride (dichloromethane), chloroform, and carbon tetrachloride are examples of common halogenated organic solvents. Dispose of all halogenated solvents in the container designated for them.

Strong inorganic acids and bases. Strong acids such as hydrochloric, sulfuric, and nitric acid will be collected in specially marked containers. Strong bases such as sodium hydroxide and potassium hydroxide will also be collected in specially designated containers.

Aqueous solutions. Aqueous solutions will be collected in a specially marked waste container. It is not necessary to separate each type of aqueous solution (unless the solution contains heavy metals); rather, unless otherwise instructed, you may combine all aqueous solutions into the same waste container. Although many types of solutions (aqueous sodium bicarbonate, aqueous sodium chloride, and so on) may seem innocuous and it may seem that their disposal down the sink drain is not likely to cause harm, many communities are becoming increasingly restrictive about what substances they will permit to enter municipal sewage-treatment systems. In light of this trend toward greater caution, it is important to develop good laboratory habits regarding the disposal of *all* chemicals.

Heavy metals. Many heavy metal ions such as mercury and chromium are highly toxic and should be disposed of in specifically designated waste containers.

Whichever method is used, the waste containers must eventually be labeled with a complete list of each substance that is present in the waste. Individual waste containers are collected, and their contents are consolidated and placed into drums for transport to the waste-disposal site. Even these drums must bear labels that detail each of the substances contained in the waste.

In either waste-handling method, certain principles will always apply:

- Aqueous solutions should not be mixed with organic liquids.
- Concentrated acids should be stored in separate containers; certainly they must *never* be allowed to come into contact with organic waste.
- Organic materials that contain halogen atoms (fluorine, chlorine, bromine, or iodine) should be stored in separate containers from those used to store materials that do not contain halogen atoms.

In each experiment in this textbook, we have suggested a method of collecting and storing wastes. Your instructor may opt to use another method for collecting wastes.

E. Use of Flames

Even though organic solvents are frequently flammable (for example, hexane, diethyl ether, methanol, acetone, and petroleum ether), there are certain laboratory procedures for which a flame must be used. Most often, these procedures involve

an aqueous solution. In fact, as a general rule, use a flame to heat only aqueous solutions. Heating methods that do not use a flame are discussed in detail in Technique 6. Most organic solvents boil below 100°C, and an aluminum block, heating mantle, sand bath, or water bath may be used to heat these solvents safely. Common organic solvents are listed in Technique 10, Table 10.3. Solvents marked in the table with boldface type will burn. Diethyl ether, pentane, and hexane are especially dangerous, because in combination with the correct amount of air, they may explode.

Some common-sense rules apply to using a flame in the presence of flammable solvents. Again, we stress that you should check to see whether anyone in your vicinity is using flammable solvents before you ignite any open flame. If someone is using a flammable solvent, move to a safer location before you light your flame. Your laboratory should have an area set aside for using a burner to prepare micro-pipettes or other pieces of glassware.

The drainage troughs or sinks should never be used to dispose of flammable organic solvents. They will vaporize if they are low boiling and may encounter a flame farther down the bench on their way to the sink.

F. Inadvertently Mixed Chemicals

To avoid unnecessary hazards of fire and explosion, never pour any reagent back into a stock bottle. There is always the chance that you may accidentally pour back some foreign substance that will react explosively with the chemical in the stock bottle. Of course, by pouring reagents back into the stock bottles, you may also introduce impurities that could spoil the experiment for the person using the stock reagent after you. Pouring reagents back into bottles is not only a dangerous practice, but also an inconsiderate one. Thus, you should not take more chemicals than you need.

G. Unauthorized Experiments

Never undertake any unauthorized experiments. The risk of an accident is high, particularly if the experiment has not been completely checked to reduce hazards. Never work alone in the laboratory. The laboratory instructor or supervisor must always be present.

H. Food in the Laboratory

Because all chemicals are potentially toxic, avoid accidentally ingesting any toxic substance; therefore, never eat or drink any food while in the laboratory. There is always the possibility that whatever you are eating or drinking may become contaminated with a potentially hazardous material.

I. Clothing

Always wear closed shoes in the laboratory; open-toed shoes or sandals offer inadequate protection against spilled chemicals or broken glass. Do not wear your best clothing in the laboratory because some chemicals can make holes in or permanent stains on your clothing. To protect yourself and your clothing, it is advisable to wear a full-length laboratory apron or coat.

When working with chemicals that are very toxic, wear some type of gloves. Disposable gloves are inexpensive, offer good protection, provide acceptable “feel,” and can be bought in many departmental stockrooms and college bookstores. Disposable latex surgical or polyethylene gloves are the least expensive type of glove; they are satisfactory when working with inorganic reagents and solutions. Better protection is afforded by disposable nitrile gloves. This type of glove

provides good protection against organic chemicals and solvents. Heavier nitrile gloves are also available.

Finally, hair that is shoulder length or longer should be tied back. This precaution is especially important if you are working with a burner.

J. First Aid: Cuts, Minor Burns, and Acid or Base Burns

If any chemical enters your eyes, immediately irrigate the eyes with copious quantities of water. Tempered (slightly warm) water, if available, is preferable. Be sure that the eyelids are kept open. Continue flushing the eyes in this way for 15 minutes.

In case of a cut, wash the wound well with water unless you are specifically instructed to do otherwise. If necessary, apply pressure to the wound to stop the flow of blood.

Minor burns caused by flames or contact with hot objects may be soothed by immediately immersing the burned area in cold water or cracked ice until you no longer feel a burning sensation. Applying salves to burns is discouraged. Severe burns must be examined and treated by a physician. For chemical acid or base burns, rinse the burned area with copious quantities of water for at least 15 minutes.

If you accidentally ingest a chemical, call the local poison control center for instructions. Do not drink anything until you have been told to do so. It is important that the examining physician be informed of the exact nature of the substance ingested.

1.2 Right-to-Know Laws

The federal government and most state governments now require that employers provide their employees with complete information about hazards in the workplace. These regulations are often referred to as **Right-to-Know Laws**. At the federal level, the Occupational Safety and Health Administration (OSHA) is charged with enforcing these regulations.

In 1990, the federal government extended the Hazard Communication Act, which established the Right-to-Know Laws, to include a provision that requires the establishment of a Chemical Hygiene Plan at all academic laboratories. Every college and university chemistry department should have a Chemical Hygiene Plan. Having this plan means that all of the safety regulations and laboratory safety procedures should be written in a manual. The plan also provides for the training of all employees in laboratory safety. Your laboratory instructor and assistants should have this training.

One of the components of Right-to-Know Laws is that employees and students have access to information about the hazards of any chemicals with which they are working. Your instructor will alert you to dangers to which you need to pay particular attention. However, you may want to seek additional information. Two excellent sources of information are labels on the bottles that come from a chemical manufacturer and **Material Safety Data Sheets (MSDSs)** or **Safety Data Sheets (SDSs)**. The MSDSs are also provided by the manufacturer and must be kept available for all chemicals used at educational institutions.

A. Material Safety Data Sheets

Reading an MSDS for a chemical can be a daunting experience, even for an experienced chemist. MSDSs contain a wealth of information, some of which must be decoded to understand. The MSDS for methanol is shown below. Only the information that might be of interest to you is described in the paragraphs that follow.

Section 1. The first part of Section 1 identifies the substance by name, formula, and various numbers and codes. Most organic compounds have more than one name. In this case, the systematic (or International Union of Pure and Applied Chemistry [IUPAC]) name is methanol, and the other names are common names or are from an older system of nomenclature. The Chemical Abstract Service Number (CAS No.) is often used to identify a substance, and it may be used to access extensive information about a substance found in many computer databases or in the library.

Section 3. The Baker SAF-T-DATA System is found on all MSDSs and bottle labels for chemicals supplied by J. T. Baker, Inc. For each category listed, the number indicates the degree of hazard. The lowest number is 0 (very low hazard), and the highest number is 4 (extreme hazard). The Health category refers to damage involved when the substance is inhaled, ingested, or absorbed. Flammability indicates the tendency of a substance to burn. Reactivity refers to how reactive a substance is with air, water, or other substances. The last category, Contact, refers to how hazardous a substance is when it comes in contact with external parts of the body. Note that this rating scale is applicable only to Baker MSDSs and labels; other rating scales with different meanings are also in common use.

Section 4. This section provides helpful information for emergency and first aid procedures.

Section 6. This part of the MSDS deals with procedures for handling spills and disposal. The information could be very helpful, particularly if a large amount of a chemical was spilled. More information about disposal is also given in Section 13.

Section 8. Much valuable information is found in Section 8. To help you understand this material, some of the more important terms used in this section are defined:

Threshold Limit Value (TLV). The American Conference of Governmental Industrial Hygienists (ACGIH) developed the TLV: This is the maximum concentration of a substance in air that a person should be exposed to on a regular basis. It is usually expressed in ppm or mg/m^3 . Note that this value assumes that a person is exposed to the substance 40 hours per week, on a long-term basis. This value may not be particularly applicable in the case of a student performing an experiment in a single laboratory period.

Permissible Exposure Limit (PEL). This has the same meaning as TLV; however, PELs were developed by OSHA. Note that for methanol, the TLV and PEL are both 200 ppm.

Section 10. The information contained in Section 10 refers to the stability of the compound and the hazards associated with mixing of chemicals. It is important to consider this information before carrying out an experiment not previously done.

Section 11. More information about the toxicity is given in this section. Another important term must first be defined:

Lethal Dose, 50% Mortality (LD_{50}). This is the dose of a substance that will kill 50% of the animals administered a single dose. Different means of administration are used, such as oral, intraperitoneal (injected into the lining of the abdominal cavity), subcutaneous (injected under the skin), and application to the surface of the skin. The LD_{50} is usually expressed in milligrams (mg) of substance per kilogram (kg) of animal weight. The lower the value of LD_{50} , the more toxic the substance. It is assumed that the toxicity in humans will be similar.

MSDS Number: M2015

Effective Date: 12/8/96

MSDS Material Safety Data Sheet

From: Mallinckrodt Baker, Inc.
222 Red School Lane
Phillipsburg, NJ 08865

MALLINCKRODT



24 Hour Emergency Telephone: 908-859-2151
CHEMTREC: 1-800-424-9300

National Response in Canada
CANUTEC: 613-996-6666

Outside U.S. and Canada
Chemtrec: 202-483-7616

NOTE: CHEMTREC, CANUTEC and National Response Center emergency numbers to be used only in the event of chemical emergencies involving a spill, leak, fire, exposure or accident involving chemicals.

All non-emergency questions should be directed to Customer Service (1-800-582-2537) for assistance.

METHYL ALCOHOL

1. Product Identification

Synonyms: Wood alcohol; methanol; carbinol
CAS No: 67-56-1
Molecular Weight: 32.04
Chemical Formula: CH₃OH
Product Codes: **J.T. Baker:**
5217, 5370, 5794, 5807, 5811, 5842, 5869, 9049, 9063, 9066, 9067, 9069, 9070, 9071, 9073, 9075, 9076, 9077, 9091, 9093, 9096, 9097, 9098, 9263, 9893
Mallinckrodt:
3004, 3006, 3016, 3017, 3018, 3024, 3041, 3701, 4295, 5160, 8814, H080, H488, H603, V079, V571

2. Composition/Information on Ingredients

Ingredient	CAS No.	Percent	Hazardous
Methyl Alcohol	67-56-1	100%	Yes

3. Hazards Identification

Emergency Overview

POISON! DANGER! VAPOR HARMFUL. MAY BE FATAL OR CAUSE BLINDNESS IF SWALLOWED. HARMFUL IF INHALED OR ABSORBED THROUGH SKIN. CANNOT BE MADE NONPOISONOUS. FLAMMABLE LIQUID AND VAPOR. CAUSES IRRITATION TO SKIN, EYES AND RESPIRATORY TRACT. AFFECTS THE LIVER.

J.T. Baker SAF-T-DATA(tm) Ratings

(Provided here for your convenience)

Health:	Flammability:	Reactivity:	Contact:
3 - Severe (Poison)	4 - Extreme (Flammable)	1 - Slight	1 - Slight
Lab Protection Equip:	GOGGLES & SHIELD; LAB COAT & APRON; VENT HOOD; PROPER GLOVES; CLASS B EXTINGUISHER		
Storage Color Code:	Red (Flammable)		

Potential Health Effects

Inhalation:

A slight irritant to the mucous membranes. Toxic effects exerted upon nervous system, particularly the optic nerve. Once absorbed into the body, it is very slowly eliminated. Symptoms of overexposure may include headache, drowsiness, nausea, vomiting, blurred vision, blindness, coma, and death. A person may get better but then worse again up to 30 hours later.

Ingestion:

Toxic. Symptoms parallel inhalation. Can intoxicate and cause blindness. Usual fatal dose: 100–125 milliliters.

Skin Contact:

Methyl alcohol is a defatting agent and may cause skin to become dry and cracked. Skin absorption can occur; symptoms may parallel inhalation exposure.

Eye Contact:

Irritant. Continued exposure may cause eye lesions.

Chronic Exposure:

Marked impairment of vision and enlargement of the liver has been reported. Repeated or prolonged exposure may cause skin irritation.

Aggravation of Pre-existing Conditions:

Persons with pre-existing skin disorders or eye problems or impaired liver or kidney function may be more susceptible to the effects of the substance.

4. First Aid Measures

Inhalation:

Remove to fresh air. If not breathing, give artificial respiration. If breathing is difficult, give oxygen. Call a physician.

Ingestion:

Induce vomiting immediately as directed by medical personnel. Never give anything by mouth to an unconscious person.

Skin Contact:

Remove any contaminated clothing. Wash skin with soap or mild detergent and water for at least 15 minutes. Get medical attention if irritation develops or persists.

Eye Contact:

Immediately flush eyes with plenty of water for at least 15 minutes, lifting lower and upper eyelids occasionally. Get medical attention immediately.

5. Fire Fighting Measures

Fire:

Flash point: 12°C (54°F) CC
Autoignition temperature: 464°C (867°F)
Flammable limits in air % by volume:
l_{el}: 7.3; u_{el}: 36
Flammable.

Explosion:

Above flash point, vapor-air mixtures are explosive within flammable limits noted above. Moderate explosion hazard and dangerous fire hazard when exposed to heat, sparks or flames. Sensitive to static discharge.

Fire Extinguishing Media:

Water spray, dry chemical, alcohol foam, or carbon dioxide.

Special Information:

In the event of a fire, wear full protective clothing and NIOSH-approved self-contained breathing apparatus with full facepiece operated in the pressure demand or other positive pressure mode. Use water spray to blanket fire, cool fire exposed containers, and to flush non-ignited spills or vapors away from fire. Vapors can flow along surfaces to distant ignition source and flash back.

6. Accidental Release Measures

Ventilate area of leak or spill. Remove all sources of ignition. Wear appropriate personal protective equipment as specified in Section 8. Isolate hazard area. Keep unnecessary and unprotected personnel from entering. Contain and recover liquid when possible. Use non-sparking tools and equipment. Collect liquid in an appropriate container or absorb with an inert material (e. g., vermiculite, dry sand, earth), and place in a chemical waste container. Do not use combustible materials, such as saw dust. Do not flush to sewer! J. T. Baker SOLUSORB® solvent adsorbent is recommended for spills of this product.

7. Handling and Storage

Protect against physical damage. Store in a cool, dry well-ventilated location, away from any area where the fire hazard may be acute. Outside or detached storage is preferred. Separate from incompatibles. Containers should be bonded and grounded for transfers to avoid static sparks. Storage and use areas should be No Smoking areas. Use non-sparking type tools and equipment, including explosion proof ventilation. Containers of this material may be hazardous when empty since they retain product residues (vapors, liquid); observe all warnings and precautions listed for the product.

8. Exposure Controls/Personal Protection

Airborne Exposure Limits:

For Methyl Alcohol:
- OSHA Permissible Exposure Limit (PEL):
200 ppm (TWA)
- ACGIH Threshold Limit Value (TLV):
200 ppm (TWA), 250 ppm (STEL) skin

Ventilation System:

A system of local and/or general exhaust is recommended to keep employee exposures below the Airborne Exposure Limits. Local exhaust ventilation is generally preferred because it can control the emissions of the contaminant at its source, preventing dispersion of it into the general work area. Please refer to the ACGIH document, "Industrial Ventilation, A Manual of Recommended Practices", most recent edition, for details.

Personal Respirator (NIOSH Approved)

If the exposure limit is exceeded, wear a supplied air, full-facepiece respirator, airlined hood, or full-facepiece self-contained breathing apparatus.

Skin Protection:

Rubber or neoprene gloves and additional protection including impervious boots, apron, or coveralls, as needed in areas of unusual exposure.

Eye Protection:

Use chemical safety goggles. Maintain eye wash fountain and quick-drench facilities in work area.

9. Physical and Chemical Properties**Appearance:**

Clear, colorless liquid.

Odor:

Characteristic odor.

Solubility:

Miscible in water.

Specific Gravity:

0.8

pH:

No information found.

% Volatiles by volume @ 21°C (70°F):

100

Boiling Point:

64.5°C (147°F)

Melting Point:

-98°C (-144°F)

Vapor Density (Air=1):

1.1

Vapor Pressure (mm Hg):

97 @ 20°C (68°F)

Evaporation Rate (BuAc=1):

5.9

10. Stability and Reactivity**Stability:**

Stable under ordinary conditions of use and storage.

Hazardous Decomposition Products:

May form carbon dioxide, carbon monoxide, and formaldehyde when heated to decomposition.

Hazardous Polymerization:

Will not occur.

Incompatibilities:

Strong oxidizing agents such as nitrates, perchlorates or sulfuric acid. Will attack some forms of plastics, rubber, and coatings. May react with metallic aluminum and generate hydrogen gas.

Conditions to Avoid:

Heat, flames, ignition sources and incompatibles.

11. Toxicological Information

Methyl Alcohol (Methanol) Oral rat LD50: 5628 mg/kg; inhalation rat LC50: 64000 ppm/4H; skin rabbit LD50: 15800 mg/kg; Irritation data-standard Draize test: skin, rabbit: 20mg/24 hr. Moderate; eye, rabbit: 100 mg/24 hr. Moderate; Investigated as a mutagen, reproductive effector.

Cancer Lists

Ingredient	--NTP Carcinogen--		IARC Category
	Known	Anticipated	
Methyl Alcohol (67-56-1)	No	No	None

12. Ecological Information**Environmental Fate:**

When released into the soil, this material is expected to readily biodegrade. When released into the soil, this material is expected to leach into groundwater. When released into the soil, this material is expected to quickly evaporate. When released into the water, this material is expected to have a half-life between 1 and 10 days. When released into water, this material is expected to readily biodegrade. When released into the air, this material is expected to exist in the aerosol phase with a short half-life. When released into the air, this material is expected to be readily degraded by reaction with photochemically produced hydroxyl radicals. When released into air, this material is expected to have a half-life between 10 and 30 days. When released into the air, this material is expected to be readily removed from the atmosphere by wet deposition.

Environmental Toxicity:

This material is expected to be slightly toxic to aquatic life.

13. Disposal Considerations

Whatever cannot be saved for recovery or recycling should be handled as hazardous waste and sent to a RCRA approved incinerator or disposed in a RCRA approved waste facility. Processing, use or contamination of this product may change the waste management options. State and local disposal regulations may differ from federal disposal regulations.

Dispose of container and unused contents in accordance with federal, state and local requirements.

14. Transport Information

Domestic (Land, D.O.T.)

Proper Shipping Name: METHANOL
Hazard Class: 3
UN/NA: UN1230 **Packing Group:** II
Information reported for product/size: 350LB

International (Water, I.M.O.)

Proper Shipping Name: METHANOL
Hazard Class: 3.2, 6.1
UN/NA: UN1230 **Packing Group:** II
Information reported for product/size: 350LB

15. Regulatory Information

Chemical Inventory Status

Ingredient	TSCA	EC	Japan	Australia	Korea	---Canada---		
						DSL	NDSL	Phil.
Methyl Alcohol (67-56-1)	Yes	Yes	Yes	Yes	Yes	Yes	No	Yes

Federal, State & International Regulations

Ingredient	---SARA 302---		---SARA 313---		CERCLA	-RCRA- 261.33	-TSCA- 8(d)
	RQ	TPQ	List	Chemical Catg.			
Methyl Alcohol (67-56-1)	No	No	Yes	No	5000	U154	No

Chemical Weapons Convention: No **TSCA 12(b):** No **CDTA:** No
SARA 311/312: Acute: Yes Chronic: Yes Fire: Yes Pressure: No Reactivity: No (Pure / Liquid)
Australian Hazchem Code: 2PE **Australian Poison Schedule:** S6
WHMIS: This MSDS has been prepared according to the hazard criteria of the Controlled Products Regulations (CPR) and the MSDS contains all of the information required by the CPR.

16. Other Information

NFPA Ratings:

Health: 1 Flammability: 3 Reactivity: 0

Label Hazard Warning:

POISON! DANGER! VAPOR HARMFUL. MAY BE FATAL OR CAUSE BLINDNESS IF SWALLOWED. HARMFUL IF INHALED OR ABSORBED THROUGH SKIN. CANNOT BE MADE NONPOISONOUS. FLAMMABLE LIQUID AND VAPOR. CAUSES IRRITATION TO SKIN, EYES AND RESPIRATORY TRACT. AFFECTS THE LIVER.

Label Precautions:

Keep away from heat, sparks and flame.
 Keep container closed.
 Use only with adequate ventilation.
 Wash thoroughly after handling.
 Avoid breathing vapor.
 Avoid contact with eyes, skin and clothing.

Label First Aid:

If swallowed, induce vomiting immediately as directed by medical personnel. Never give anything by mouth to an unconscious person. In case of contact, immediately flush eyes or skin with plenty of water for at least 15 minutes while removing contaminated clothing and shoes. Wash clothing before reuse. If inhaled, remove to fresh air. If not breathing give artificial respiration. If breathing is difficult, give oxygen. In all cases get medical attention immediately.

Product Use:

Laboratory Reagent.

Revision Information:

New 16 section MSDS format, all sections have been revised.


Disclaimer:

Mallinckrodt Baker, Inc. provides the information contained herein in good faith but makes no representation as to its comprehensiveness or accuracy. This document is intended only as a guide to the appropriate precautionary handling of the material by a properly trained person using this product. Individuals receiving the information must exercise their independent judgment in determining its appropriateness for a particular purpose. MALLINCKRODT BAKER, INC. MAKES NO REPRESENTATIONS OR WARRANTIES, EITHER EXPRESS OR IMPLIED, INCLUDING WITHOUT LIMITATION ANY WARRANTIES OR MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE WITH RESPECT TO THE INFORMATION SET FORTH HEREIN OR THE PRODUCT TO WHICH THE INFORMATION REFERS. ACCORDINGLY, MALLINCKRODT BAKER, INC. WILL NOT BE RESPONSIBLE FOR DAMAGES RESULTING FROM USE OF OR RELIANCE UPON THIS INFORMATION.

Prepared By: Strategic Services Division
 Phone Number: (314) 539-1600 (U.S.A.)

Unless you have considerably more knowledge about chemical toxicity, the information in Sections 8 and 11 is most useful for comparing the toxicity of one substance with another. For example, the TLV for methanol is 200 ppm, whereas the TLV for benzene is 10 ppm. Clearly, performing an experiment involving benzene would require much more stringent precautions than an experiment involving methanol. One of the LD₅₀ values for methanol is 5628 mg/kg. The comparable LD₅₀ value of aniline is 250 mg/kg. Clearly, aniline is much more toxic, and because it is easily absorbed through the skin, it presents a significant hazard. It should also be mentioned that both TLV and PEL ratings assume that the worker comes in contact with a substance on a repeated and long-term basis. Thus, even if a chemical has a relatively low TLV or PEL, it does not mean that using it for one experiment will present a danger to you. Furthermore, by performing experiments using small amounts of chemicals and with proper safety precautions, your exposure to organic chemicals in this course will be minimal.

Section 16. Section 16 contains the National Fire Protection Association (NFPA) rating. This is similar to the Baker SAF-T-DATA (discussed in Section 3), except that the number represents the hazards when a fire is present. The order here is Health, Flammability, and Reactivity. Often, this is presented in graphic form on a label (see figure). The small diamonds are often color coded: blue for Health, red for Flammability, and yellow for Reactivity. The bottom diamond (white) is sometimes used to display graphic symbols denoting unusual reactivity, hazards, or special precautions to be taken.

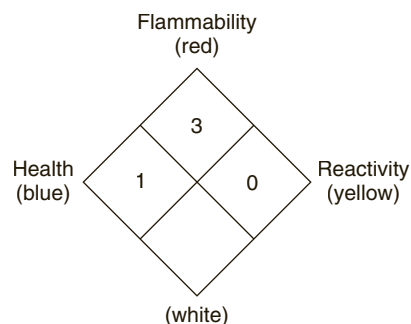

A newer version of the NFPA diamond has come into use since about 2001. This is the Hazardous Materials Identification System (HMIS) color bar, developed by the American Coatings Association, working with OSHA, to augment the older NFPA warning system. The HMIS system uses a color bar.

Chemical Name	
HEALTH	0
FLAMMABILITY	0
PHYSICAL HAZARD	0
PERSONAL PROTECTION	0

The uppermost bar, which appears in blue, represents the degree of health hazard associated with a substance. A rating of zero indicates that there is no significant health risk, while a rating of four indicates a major, life-threatening risk.

An asterisk may also appear, indicating that chronic exposure could also cause health problems.

The second, red bar rates the degree to which a substance is flammable. The third, orange bar is a hazard rating. The final bar, which is white, indicates the level of personal protection that is recommended while working with the substance. The white bar used in the HMIS system is significantly different than the white area used in the NFPA diamond.

B. Safety Data Sheets – SDS and GHS

Since the advent of the Material Safety Data Sheet, a different reporting style has evolved. In 2012, OSHA's Hazard Communications Standards (HCS) were upgraded. The purpose was to mandate a more consistent format and to align the format with the Globally Harmonized System (GHS) of Classification and labelling of Chemicals which was adopted internationally. The goal was to make the SDS format more consistent with equivalent safety sheets used in other countries. At the same time, the MSDS was renamed to a simpler SDS (Safety Data Sheet) acronym.

The format and content of the newer SDS is essentially the same as the older MSDS. In practice, there is no significant difference between the information found on an MSDS and what appears on an SDS. The SDS, however, makes greater use of pictograms to highlight some critical safety information. The following information is taken directly from the OSHA "Quick Card" that explains the new SDS requirements:

"The new OSHA Hazard Communication Standard (HCS) requires chemical manufacturers, distributors, or importers to provide Safety Data Sheets (SDSs) (formerly known as Material Safety Data Sheets or MSDSs) to communicate the hazards of hazardous chemical products. As of June 1, 2015, the HCS will require new SDSs to be in a uniform format, and include the section numbers, the headings, and associated information under the headings below:

Section 1, Identification includes product identifier; manufacturer or distributor name, address, phone number; emergency phone number; recommended use; restrictions on use.

Section 2, Hazard(s) identification includes all hazards regarding the chemical; required label elements. [This section requires use of the safety pictograms.]

Section 3, Composition/information on ingredients includes information on chemical ingredients; trade secret claims.

Section 4, First-aid measures includes important symptoms/effects, acute, delayed; required treatment.

Section 5, Fire-fighting measures lists suitable extinguishing techniques, equipment; chemical hazards from fire.

Section 6, Accidental release measures lists emergency procedures; protective equipment; proper methods of containment and cleanup.

Section 7, Handling and storage lists precautions for safe handling and storage, including incompatibilities.

Section 8, Exposure controls/personal protection lists OSHA's Permissible Exposure Limits (PELs); ACGIH Threshold Limit Values (TLVs); and any other exposure limit used or recommended by the chemical manufacturer, importer, or employer preparing the SDS where available as well as appropriate engineering controls; personal protective equipment (PPE).

Section 9, Physical and chemical properties lists the chemical's characteristics.

Section 10, Stability and reactivity lists chemical stability and possibility of hazardous reactions.

Section 11, Toxicological information includes routes of exposure; related symptoms, acute and chronic effects; numerical measures of toxicity.

Section 12, Ecological information

Section 13, Disposal considerations










Section 14, Transport information

Section 15, Regulatory information

Section 16, Other information, includes the date of preparation or last revision.

Employers must ensure that SDSs are readily accessible to employees.”

HCS Pictograms and Hazards

<p>Health Hazard</p>  <ul style="list-style-type: none"> • Carcinogen • Mutagen • Reproductive Toxicity • Aspiration Toxicity • Target Organ Toxicity • Respiratory Sensitizer 	<p>Fire Hazard</p>  <ul style="list-style-type: none"> • Flammable • Pyrophoric • Emits Flammable Gas • Self-Reacting • Self-Heating • Organic Peroxide 	<p>Irritant</p>  <ul style="list-style-type: none"> • Irritant (skin and eye) • Irritant (respiratory) • Acute Toxicity (harmful) • Skin Sensitizer • Narcotic Effects
<p>Pressurized</p>  <ul style="list-style-type: none"> • Gases Under Pressure • Explosion Hazard 	<p>Corrosive</p>  <ul style="list-style-type: none"> • Skin Corrosion/ Burns • Eye Damage • Corrosive to Metals 	<p>Explosive</p>  <ul style="list-style-type: none"> • Explosive • Self-Reacting • Organic Peroxide
<p>Oxidizer</p>  <ul style="list-style-type: none"> • Oxidizer • Causes Fires 	<p>Environmental Hazard</p>  <ul style="list-style-type: none"> • Aquatic Toxicity 	<p>Poison</p>  <ul style="list-style-type: none"> • Acute Toxicity (fatal or toxic)

C. Bottle Labels

Reading the label on a bottle can be a very helpful way of learning about the hazards of a chemical. The amount of information varies greatly, depending on which company supplied the chemical.

Apply some common sense when you read MSDSs and bottle labels. Using these chemicals does not mean you will experience the consequences that can potentially result from exposure to each chemical. For example, an MSDS for sodium chloride states, "Exposure to this product may have serious adverse health effects." Despite the apparent severity of this cautionary statement, it would not be reasonable to expect people to stop using sodium chloride in a chemistry experiment or to stop sprinkling a small amount of it (as table salt) on eggs to enhance their flavor. In many cases, the consequences described in MSDSs from exposure to chemicals are somewhat overstated, particularly for students using these chemicals to perform a laboratory experiment.

1.3 Common Solvents

Most organic chemistry experiments involve an organic solvent at some step in the procedure. A list of common organic solvents follows, with a discussion of toxicity, possible carcinogenic properties, and precautions that you should use when handling these solvents. A tabulation of the compounds currently suspected of being carcinogens appears at the end of Technique 1.

Acetic Acid. Glacial acetic acid is corrosive enough to cause serious acid burns on the skin. Its vapors can irritate the eyes and nasal passages. Care should be exercised not to breathe the vapors and not to allow them to escape into the laboratory.

Acetone. Relative to other organic solvents, acetone is not very toxic. It is flammable, however. Do not use acetone near open flames.

Benzene. Benzene can damage bone marrow, it causes various blood disorders, and its effects may lead to leukemia. Benzene is considered a serious carcinogenic hazard. It is absorbed rapidly through the skin and also poisons the liver and kidneys. In addition, benzene is flammable. Because of its toxicity and its carcinogenic properties, benzene should not be used in the laboratory; you should use some less dangerous solvent instead. Toluene is considered a safer alternative solvent in procedures that specify benzene.

Carbon Tetrachloride. Carbon tetrachloride can cause serious liver and kidney damage, as well as skin irritation and other problems. It is absorbed rapidly through the skin. In high concentrations, it can cause death as a result of respiratory failure. Moreover, carbon tetrachloride is suspected of being a carcinogenic material. Although this solvent has the advantage of being nonflammable (in the past, it was used on occasion as a fire extinguisher), it can cause health problems, so it should not be used routinely in the laboratory. If no reasonable substitute exists, however, it must be used in small quantities, as in preparing samples for infrared (IR) and nuclear magnetic resonance (NMR) spectroscopy. In such cases, you must use it in a hood.

Chloroform. Chloroform is similar to carbon tetrachloride in its toxicity. It has been used as an anesthetic. However, chloroform is currently on the list of suspected carcinogens. Because of this, do not use chloroform routinely as a solvent in the laboratory. If it is occasionally necessary to use chloroform as a solvent for special samples, then you must use it in a hood. Methylene chloride is usually found to be a safer substitute in procedures that specify chloroform as a solvent. Deuteriochloroform, CDCl_3 , is a common solvent for NMR spectroscopy. Caution dictates that you should treat it with the same respect as chloroform.

1,2-Dimethoxyethane (Ethylene Glycol Dimethyl Ether or Monoglyme). Because it is miscible with water, 1,2-dimethoxyethane is a useful alternative to solvents such as dioxane and tetrahydrofuran, which may be more hazardous. 1,2-Dimethoxyethane is flammable and should not be handled near an open flame. Upon long exposure of 1,2-dimethoxyethane to light and oxygen, explosive peroxides may form. 1,2-Dimethoxyethane is also a possible reproductive toxin.

Dioxane. Dioxane has been used widely because it is a convenient, water-miscible solvent. It is now suspected, however, of being carcinogenic. It is also toxic, affecting the central nervous system, liver, kidneys, skin, lungs, and mucous membranes. Dioxane is also flammable and tends to form explosive peroxides when it is exposed to light and air. Because of its carcinogenic properties, it is no longer used in the laboratory unless absolutely necessary. Either 1,2-dimethoxyethane or tetrahydrofuran is a suitable, water-miscible alternative solvent.

Ethanol. Ethanol has well-known properties as an intoxicant. In the laboratory, the principal danger arises from fires, because ethanol is a flammable solvent. When using ethanol, take care to work where there are no open flames.

Ether (diethyl ether). The principal hazard associated with diethyl ether is fire or explosion. Ether is probably the most flammable solvent found in the laboratory. Because ether vapors are much denser than air, they may travel along a laboratory bench for a considerable distance from their source before being ignited. Before using ether, it is very important to be sure that no one is working with matches or any open flame. Ether is not a particularly toxic solvent, although in high enough concentrations it can cause drowsiness and perhaps nausea. It has been used as a general anesthetic. Ether can form highly explosive peroxides when exposed to air. Consequently, you should never distill it to dryness.

Hexane. Hexane may be irritating to the respiratory tract. It can also act as an intoxicant and a depressant of the central nervous system. It can cause skin irritation because it is an excellent solvent for skin oils. The most serious hazard, however, comes from its flammability. The precautions recommended for using diethyl ether in the presence of open flames apply equally to hexane.

Ligroin. See Hexane.

Methanol. Much of the material outlining the hazards of ethanol applies to methanol. Methanol is more toxic than ethanol; ingestion can cause blindness and even death. Because methanol is more volatile, the danger of fires is more acute.

Methylene Chloride (Dichloromethane). Methylene chloride is not flammable. Unlike other members of the class of chlorocarbons, it is not currently considered a serious carcinogenic hazard. Recently, however, it has been the subject of much serious investigation, and there have been proposals to regulate it in industrial situations in which workers have high levels of exposure on a day-to-day basis. Methylene chloride is less toxic than chloroform and carbon tetrachloride. It can cause liver damage when ingested, however, and its vapors may cause drowsiness or nausea.

Pentane. See Hexane.

Petroleum Ether. See Hexane.

Pyridine. Some fire hazard is associated with pyridine. However, the most serious hazard arises from its toxicity. Pyridine may depress the central nervous system;

irritate the skin and respiratory tract; damage the liver, kidneys, and gastrointestinal system; and even cause temporary sterility. You should treat pyridine as a highly toxic solvent and handle it only in the fume hood.

Tetrahydrofuran. Tetrahydrofuran may cause irritation of the skin, eyes, and respiratory tract. It should never be distilled to dryness because it tends to form potentially explosive peroxides on exposure to air. Tetrahydrofuran does present a fire hazard.

Toluene. Unlike benzene, toluene is not considered a carcinogen. However, it is at least as toxic as benzene. It can act as an anesthetic and damage the central nervous system. If benzene is present as an impurity in toluene, expect the usual hazards associated with benzene. Toluene is also a flammable solvent, and the usual precautions about working near open flames should be applied.

You should not use certain solvents in the laboratory because of their carcinogenic properties. Benzene, carbon tetrachloride, chloroform, and dioxane are among these solvents. For certain applications, however, notably as solvents for infrared or NMR spectroscopy, there may be no suitable alternative. When it is necessary to use one of these solvents, use safety precautions and refer to the discussions in Techniques 25–28.

Because relatively large amounts of solvents may be used in a large organic laboratory class, your laboratory supervisor must take care to store these substances safely. Only the amount of solvent needed for a particular experiment should be kept in the laboratory. The preferred location for bottles of solvents being used during a class period is in a hood. When the solvents are not being used, they should be stored in a fireproof storage cabinet for solvents. If possible, this cabinet should be ventilated into the fume hood system.

1.4 Carcinogenic Substances

A **carcinogen** is a substance that causes cancer in living tissue. The usual procedure for determining whether a substance is carcinogenic is to expose laboratory animals to high dosages over a long period. It is not clear whether short-term exposure to these chemicals carries a comparable risk, but it is prudent to use these substances with special precautions.

Many regulatory agencies have compiled lists of carcinogenic substances or substances suspected of being carcinogenic. Because these lists are inconsistent, compiling a definitive list of carcinogenic substances is difficult. The following common substances are included in many of these lists.

Acetamide	4-Methyl-2-oxetanone (β -butyrolactone)
Acrylonitrile	1-Naphthylamine
Asbestos	2-Naphthylamine
Benzene	N-Nitroso compounds
Benzidine	2-Oxetanone (β -propiolactone)
Carbon tetrachloride	Phenacetin
Chloroform	Phenylhydrazine and its salts
Chromic oxide	Polychlorinated biphenyl (PCB)
Coumarin	Progesterone
Diazomethane	Styrene oxide
1,2-Dibromoethane	Tannins

Dimethyl sulfate	Testosterone
<i>p</i> -Dioxane	Thioacetamide
Ethylene oxide	Thiourea
Formaldehyde	<i>o</i> -Toluidine
Hydrazine and its salts	Trichloroethylene
Lead (II) acetate	Vinyl chloride

REFERENCES

- Aldrich Catalog and Handbook of Fine Chemicals*; Aldrich Chemical: Milwaukee, WI, current edition.
- Armour, M. A. *Pollution Prevention and Waste Minimization in Laboratories*; Reinhardt, P. A., Leonard, K. L., Ashbrook, P. C., Eds.; Lewis Publishers: Boca Raton, FL, 1996.
- Colonna, G. R., Ed. *Fire Protection Guide on Hazardous Materials*; National Fire Protection Association: Quincy, MA, 2010.
- Flinn Chemical Catalog Reference Manual*; Flinn Scientific: Batavia, IL, 2014.
- Gosselin, R. E.; Smith, R. P.; Hodge, H. C. *Clinical Toxicology of Commercial Products*, 5th ed.; Williams & Wilkins: Baltimore, MD, 1984.
- Identifying and Evaluating Hazards in Research Laboratories*; Committee on Chemical Safety, American Chemical Society: Washington, DC, 2015.
- Lenga, R. E., Ed. *The Sigma-Aldrich Library of Chemical Safety Data*, 2nd ed.; Sigma-Aldrich: Milwaukee, WI, 1988.
- Lewis, R. J. *Carcinogenically Active Chemicals: A Reference Guide*; Van Nostrand Reinhold: New York, 1990.
- Lewis, R. J. *Sax's Dangerous Properties of Industrial Materials*, 11th ed; Van Nostrand Reinhold: New York, 2007.
- The Merck Index*, 15th ed; Merck and Co.: Rahway, NJ, 2013.
- Prudent Practices in the Laboratory: Handling and Disposal of Chemicals*; Committee on Prudent Practices for Handling, Storage, and Disposal of Chemicals in Laboratories; Board on Chemical Sciences and Technology; Commission on Physical Sciences, Mathematics, and Applications; National Research Council, National Academy Press: Washington, DC, 2011.
- Renfrew, M. M., Ed. *Safety in the Chemical Laboratory*; Division of Chemical Education, American Chemical Society: Easton, PA, 1967–1991.
- Safe Science: Promoting a Culture of Safety in Academic Chemical Research*; National Research Council, National Academy Press: Washington, DC, 2014.
- Safety in Academic Chemistry Laboratories*, 7th ed.; Committee on Chemical Safety, American Chemical Society: Washington, DC, 2003.
- Sax, N. I., Lewis, R. J., Eds. *Rapid Guide to Hazardous Chemicals in the Work Place*, 4th ed.; Van Nostrand Reinhold: New York, 2000.

Useful Safety-Related Internet Addresses

Interactive Learning Paradigms, Inc.

<http://www.ilpi.com/msds/>

This is an excellent general site for MSDS sheets. The site lists chemical manufacturers and suppliers. Selecting a company will take you directly to the appropriate place to obtain an MSDS sheet. Many of the sites listed require you to register in order to obtain an MSDS sheet for a particular chemical. Ask your departmental or college safety supervisor to obtain the information for you.

Acros chemicals and Fisher Scientific

<https://www1.fishersci.com/>

Alfa Aesar

<http://www.alfa.com/en/advanced-search>

Cornell University, Department of Environmental Health and Safety

<https://sp.ehs.cornell.edu/lab-research-safety/research-safety/msds/Pages/default.aspx>

This is an excellent searchable database of more than 325,000 MSDS files. No registration is required.

Eastman Kodak

http://www.kodak.com/ek/US/en/About_Kodak_Top/Global_Sustainability/Material_Safety_Data_Sheets.htm

<http://www.emdmillipore.com/US/en/support/safety/safety-data-sheets/Ivmb.qB.TzsAAAFcXd4Xr74u.nav>

http://www.emdchemicals.com/corporate/emd_corporate.asp

J. T. Baker and Mallinckrodt Laboratory Chemicals

<http://www.jtbaker.com/asp/Catalog.asp>

The Merck Index Online has many useful features including structure searches. <https://www.rsc.org/merck-index>

National Institute for Occupational Safety and Health (NIOSH) has an excellent "Website" that includes databases and information resources, including links:

<http://www.cdc.gov/niosh>

Sigma, Aldrich and Fluka

<http://www.sigmaaldrich.com/safety-center.html>

VWR Scientific Products

<https://us.vwr.com/store/search/searchMSDS.jsp>

The Laboratory Notebook, Calculations, and Laboratory Records

In the Introduction to this book, we mentioned the importance of advance preparation for laboratory work. Presented here are some suggestions about what specific information you should try to obtain in your advance studying. Because much of this information must be obtained while preparing your laboratory notebook, the two subjects, advance study and notebook preparation, are developed simultaneously.

An important part of any laboratory experience is learning to maintain very complete records of every experiment undertaken and every item of data obtained. Far too often, careless recording of data and observations has resulted in mistakes, frustration, and lost time due to needless repetition of experiments. If reports are required, you will find that proper collection and recording of data can make your report writing much easier.

Because organic reactions are seldom quantitative, special problems result. Frequently, reagents must be used in large excess to increase the amount of product. Some reagents are expensive, and, therefore, care must be used in measuring the amounts of these substances. Very often, many more reactions take place than you desire. These extra reactions, or **side reactions**, may form products other than the desired product. These are called **side products**. For all of these reasons, you must plan your experimental procedure carefully before undertaking the actual experiment.

2.1 The Notebook

For recording data and observations during experiments, use a *bound notebook*. The notebook should have consecutively numbered pages. If it does not, number the pages immediately. A spiral-bound notebook or any other notebook from which the pages can be removed easily is not acceptable, because the possibility of losing the pages is great.

All data and observations must be recorded in the notebook. Paper towels, napkins, toilet tissue, or scratch paper tend to become lost or destroyed. It is bad laboratory practice to record information on such random and perishable pieces of paper. All entries must be recorded in *permanent ink*. It can be frustrating to have important information disappear from the notebook because it was recorded in washable ink or pencil and could not survive a flood caused by the student at the next position on the bench. Because you will be using your notebook in the laboratory, the book will probably become soiled or stained by chemicals, filled with scratched-out entries, or even slightly burned. That is expected and is a normal part of laboratory work.



Your instructor may check your notebook at any time, so you should always have it up to date. If your instructor requires reports, you can prepare them quickly from the material recorded in the laboratory notebook.

2.2 Notebook Format A. Advance Preparation

Individual instructors vary greatly in the type of notebook format they prefer; such variation stems from differences in philosophies and experience. You must obtain specific directions from your own instructor for preparing a notebook. Certain features, however, are common to most notebook formats. The following discussion indicates what might be included in a typical notebook.

It will be very helpful and you can save much time in the laboratory if for each experiment you know the main reactions, the potential side reactions, the mechanism, and the stoichiometry, and you understand fully the procedure and the theory underlying it before you come to the laboratory. Understanding the procedure by which the desired product is to be separated from undesired materials is also very important. If you examine each of these topics before coming to class, you will be prepared to do the experiment efficiently. You will have your equipment and reagents already prepared when they are to be used. Your reference material will be at hand when you need it. Finally, with your time efficiently organized, you will be able to take advantage of long reaction or reflux periods to perform other tasks, such as doing shorter experiments or finishing previous ones.


For experiments in which a compound is synthesized from other reagents, that is, **preparative experiments**, it is essential to know the main reaction. To perform stoichiometric calculations, you should balance the equation for the main reaction. Therefore, before you begin the experiment, your notebook should contain the balanced equation for the pertinent reaction. Using the preparation of isopentyl acetate, or banana oil, as an example, you should write the following:

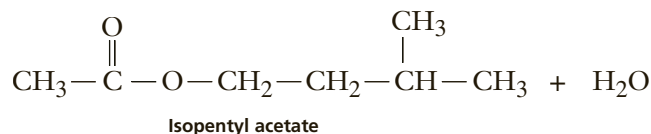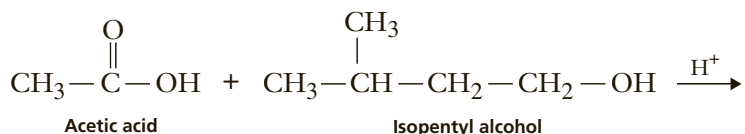

Also, before beginning the experiment enter in the notebook the possible side reactions that divert reagents into contaminants (side products). You will have to separate these side products from the major product during purification.

You should list physical constants such as melting points, boiling points, densities, and molecular weights in the notebook when this information is needed to perform an experiment or to do calculations. These data are located in sources such as the *CRC Handbook of Chemistry and Physics*, *The Merck Index*, *Lange's Handbook of Chemistry*, or the *Aldrich Handbook of Fine Chemicals*. Write physical constants required for an experiment in your notebook before you come to class.

Advance preparation may also include examining some subjects, information not necessarily recorded in the notebook, that should prove useful in understanding the experiment. Included among these subjects are an understanding of the mechanism of the reaction, an examination of other methods by which the same compound might be prepared, and a detailed study of the experimental procedure. Many students find that an outline of the procedure, prepared *before* they come to class, helps them use their time more efficiently once they begin the experiment. Such an outline could very well be prepared on some loose sheet of paper rather than in the notebook itself.

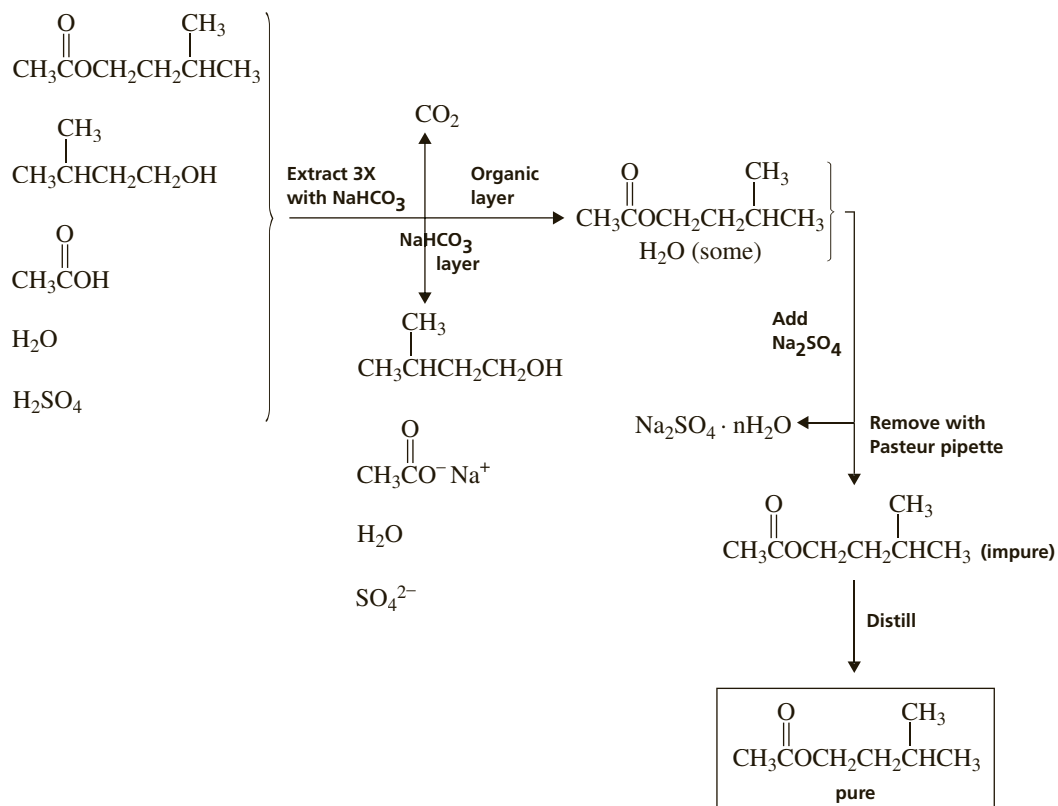**Figure 2.1**

Separation scheme for isopentyl acetate.

Once the reaction has been completed, the desired product does not magically appear as purified material; it must be isolated from a frequently complex mixture of side products, unreacted starting materials, solvents, and catalysts. You should try to outline a **separation scheme** in your notebook for isolating the product from its contaminants. At each stage, you should try to understand the reason for the particular instruction given in the experimental procedure. This not only will familiarize you with the basic separation and purification techniques used in organic chemistry but also will help you understand when to use these techniques. Such an outline might take the form of a flowchart. For example, see the separation scheme for isopentyl acetate (see Figure 2.1). Careful attention to understanding the separation, besides familiarizing you with the procedure by which the desired product is separated from impurities in your particular experiments, may prepare you for original research in which no experimental procedure exists.

In designing a separation scheme, note that the scheme outlines those steps undertaken once the reaction period has been concluded. For this reason, the represented scheme does not include steps such as the addition of the reactants (isopentyl alcohol and acetic acid) and the catalyst (sulfuric acid) or the heating of the reaction mixture. See Technique 12, Section 12.12, for a more thorough discussion of how to outline a separation scheme using a flowchart.

For experiments in which a compound is isolated from a particular source and is not prepared from other reagents, some information described in this section will not be applicable. Such experiments are called **isolation experiments**. A typical isolation experiment involves isolating a pure compound from a natural source. Examples include isolating caffeine from tea or isolating cinnamaldehyde from

cinnamon. Although isolation experiments require somewhat different advance preparation, this advance study may include looking up physical constants for the compound isolated and outlining the isolation procedure. A detailed examination of the separation scheme is very important here because it is the heart of such an experiment.

B. Laboratory Records

When you begin the actual experiment, keep your notebook nearby so you will be able to record those operations you perform. When working in the laboratory, your notebook serves as a place in which to record a rough transcript of your experimental method. Data from actual weighings, volume measurements, and determinations of physical constants are also noted. This section of your notebook should *not* be prepared in advance. The purpose is not to write a recipe but rather to record what you *did* and what you *observed*. These observations will help you write reports without resorting to memory. They will also help you or other workers repeat the experiment in as nearly as possible the same way. The sample notebook pages found in Figures 2.2 and 2.3 illustrate the type of data and observations that should be written in your notebook.

When your product has been prepared and purified, or isolated if it is an isolation experiment, record pertinent data such as the melting point or boiling point of the substance, its density, its index of refraction, and the conditions under which spectra were determined.

C. Calculations


A chemical equation for the overall conversion of the starting materials to products is written on the assumption of simple ideal stoichiometry. Actually, this assumption is seldom realized. Side reactions or competing reactions will also occur, giving other products. For some synthetic reactions, an equilibrium state will be reached in which an appreciable amount of starting material is still present and can be recovered. Some of the reactant may also remain if it is present in excess or if the reaction was incomplete. A reaction involving an expensive reagent illustrates another reason for needing to know how far a particular type of reaction converts reactants to products. In such a case, it is preferable to use the most efficient method for this conversion. Thus, information about the efficiency of conversion for various reactions is of interest to the person contemplating the use of these reactions.

The quantitative expression for the efficiency of a reaction is found by calculating the **yield** for the reaction. The **theoretical yield** is the number of grams of the product expected from the reaction on the basis of ideal stoichiometry, with side reactions, reversibility, and losses ignored. To calculate the theoretical yield, it is first necessary to determine the **limiting reagent**. The limiting reagent is the reagent that is not present in excess and on which the overall yield of product depends. The method for determining the limiting reagent in the isopentyl acetate experiment is illustrated in the sample notebook pages shown in Figures 2.2 and 2.3. You should consult your general chemistry textbook for more complicated examples. The theoretical yield is then calculated from the expression:

Theoretical yield = (moles of limiting reagent)(ratio)(molecular weight of product)

The ratio here is the stoichiometric ratio of product to limiting reagent. In preparing isopentyl acetate, that ratio is 1:1. One mole of isopentyl alcohol, under ideal circumstances, should yield 1 mole of isopentyl acetate.

THE PREPARATION OF ISOPENTYLACETATE (BANANA OIL)

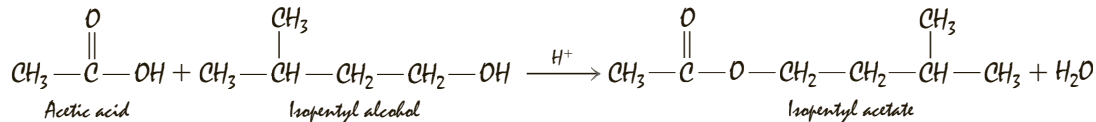Main ReactionTable of Physical Constants

	MW	BP	Density
Isopentyl alcohol	88.2	132.8°C	0.813 g/ml
Acetic acid	60.1	118	1.06
Isopentyl acetate	130.2	142	0.876

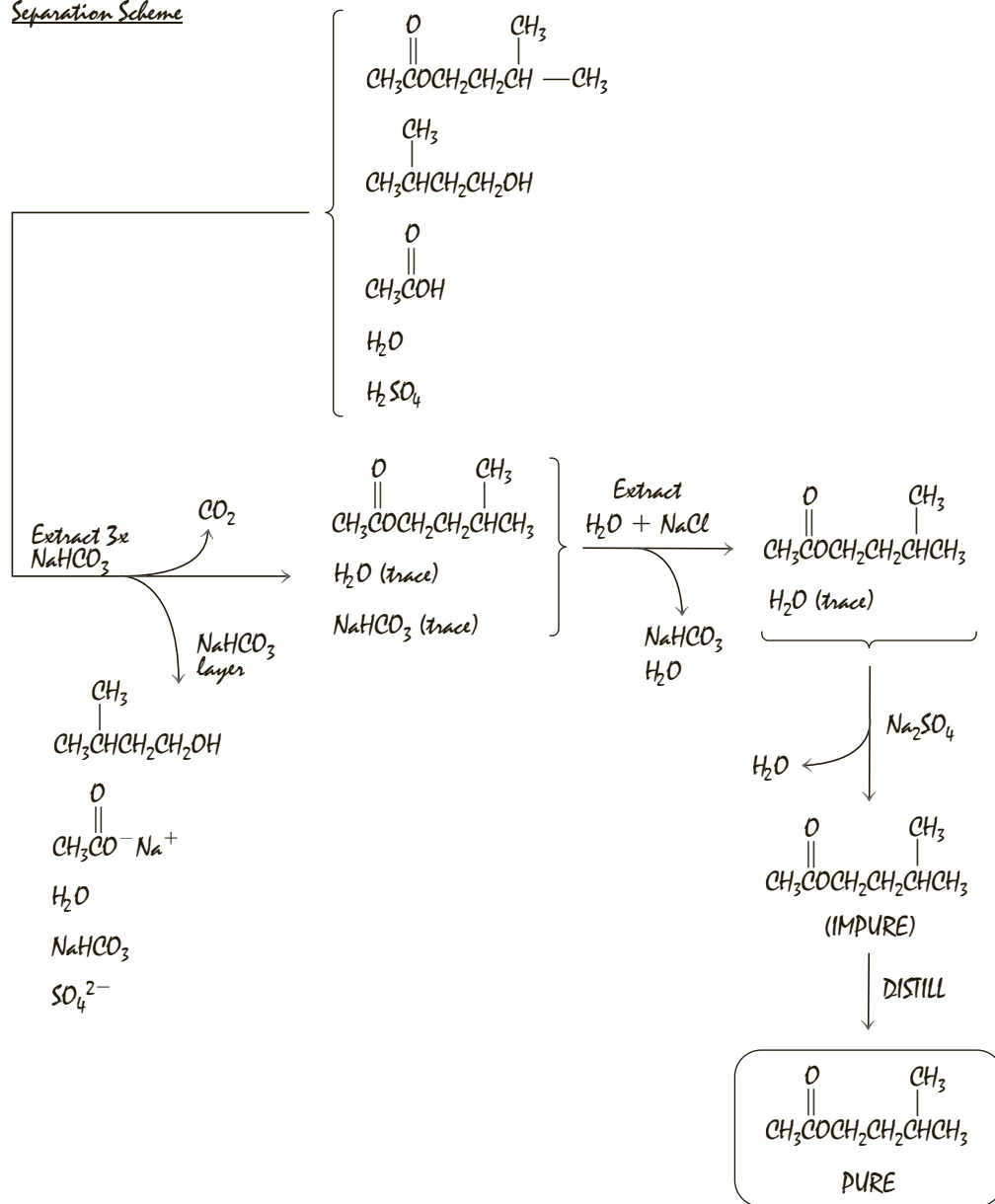
Separation Scheme

Figure 2.2

A sample notebook, page 1.

Data and Observations

7.5 mL of isopentyl alcohol was added to a pre-weighed 50-mL round-bottomed flask:

Flask + alcohol	139.75 g
Flask	133.63 g
	6.12 g isopentyl alcohol

Glacial acetic acid (10 mL) and 2 mL of concentrated sulfuric acid were also added to the flask, with swirling, along with several boiling stones. A water-cooled condenser was attached to the flask. The reaction was allowed to boil, using a heating mantle, for about one hour. The color of the reaction mixture was brownish-yellow.

After the reaction mixture had cooled to room temperature, the boiling stones were removed, and the reaction mixture was poured into a separatory funnel. About 30 mL of cold water was added to the separatory funnel. The reaction flask was rinsed with 5 mL of cold water, and the water was also added to the separatory funnel. The separatory funnel was shaken, and the lower aqueous layer was removed and discarded. The organic layer was extracted twice with two 10–15-mL portions of 5% aqueous sodium bicarbonate. During the first extraction, much CO₂ was given off, but the amount of gas evolved was markedly diminished during the second extraction. The organic layer was a light yellow in color. After the second extraction, the aqueous layer turned red litmus blue. The bicarbonate layers were discarded, and the organic layer was extracted with a 10–15-mL portion of water. A 2–3 mL portion of saturated sodium chloride solution was added during this extraction. When the aqueous layer had been removed, the upper, organic phase was transferred to a 15-mL Erlenmeyer flask. 2 g of anhydrous magnesium sulfate was added. The flask was stoppered, swirled gently, and allowed to stand for 15 mins.

The product was transferred to a 25-mL round-bottomed flask, and it was distilled by simple distillation. The distillation continued until no liquid could be observed dripping into the collection flask. After the distillation, the ester was transferred to a pre-weighed sample vial.

Sample vial + product	9.92 g
Sample vial	6.11 g
	3.81 g isopentyl acetate

The product was colorless and clear. The observed boiling point obtained during the distillation, was 140°C. An IR spectrum was obtained of the product.

Calculations

Determine limiting reagent:

$$\text{isopentyl alcohol } 6.12 \text{ g} \left(\frac{1 \text{ mol isopentyl alcohol}}{88.2 \text{ g}} \right) = 6.94 \times 10^{-2} \text{ mol}$$

$$\text{acetic acid: } (10 \text{ mL}) \left(\frac{1.06 \text{ g}}{\text{mL}} \right) \left(\frac{1 \text{ mol acetic acid}}{60.1 \text{ g}} \right) = 1.76 \times 10^{-1} \text{ mol}$$

Since they react in a 1:1 ratio, isopentyl alcohol is the limiting reagent. Theoretical yield:

$$\begin{aligned} (6.94 \times 10^{-2} \text{ mol isopentyl alcohol}) & \left(\frac{1 \text{ mol isopentyl}}{1 \text{ mol isopentyl alcohol}} \right) \left(\frac{130.2 \text{ g isopentyl acetate}}{1 \text{ mol isopentyl acetate}} \right) \\ & = 9.03 \text{ g isopentyl acetate} \end{aligned}$$

$$\text{Percentage yield} = \frac{3.81 \text{ g}}{9.03 \text{ g}} \times 100 = 42.2\%$$

Figure 2.3

A sample notebook, page 2.

The **actual yield** is simply the number of grams of desired product obtained. The **percentage yield** describes the efficiency of the reaction and is determined by

$$\text{Percentage yield} = \frac{\text{Actual yield}}{\text{Theoretical yield}} \times 100$$

Calculation of the theoretical yield and percentage yield can be illustrated using hypothetical data for the isopentyl acetate preparation:

$$\begin{aligned} \text{Theoretical yield} &= (6.94 \times 10^{-2} \text{ mol isopentyl alcohol}) \left(\frac{1 \text{ mol isopentyl acetate}}{1 \text{ mol isopentyl alcohol}} \right) \\ &\times \left(\frac{130.2 \text{ g isopentyl acetate}}{1 \text{ mol isopentyl acetate}} \right) = 9.03 \text{ g isopentyl acetate} \end{aligned}$$

$$\text{Actual yield} = 3.81 \text{ g isopentyl acetate}$$

$$\text{Percentage yield} = \frac{3.81 \text{ g}}{9.03 \text{ g}} \times 100 = 42.2\%$$

For experiments that have the principal objective of isolating a substance such as a natural product rather than preparing and purifying some reaction product, the **weight percentage recovery** and not the percentage yield is calculated. This value is determined by

$$\text{Weight percentage recovery} = \frac{\text{Weight of substance isolated}}{\text{Weight of original material}} \times 100$$

Thus, for instance, if 0.014 g of caffeine was obtained from 2.3 g of tea, the weight percentage recovery of caffeine would be

$$\text{Weight percentage recovery} = \frac{0.014 \text{ g caffeine}}{2.3 \text{ g tea}} \times 100 = 0.61\%$$

2.3 Laboratory Reports

Various formats for reporting the results of the laboratory experiments may be used. You may write the report directly in your notebook in a format similar to the sample notebook pages included in this section. Alternatively, your instructor may require a more formal report that is not written in your notebook. When you do original research, these reports should include a detailed description of all the experimental steps undertaken. Frequently, the style used in scientific periodicals such as *Journal of the American Chemical Society* is applied to writing laboratory reports. Your instructor is likely to have his or her own requirements for laboratory reports and should describe the requirements to you.

2.4 Submission of Samples

In all preparative experiments and in some isolation experiments, you will be required to submit to your instructor the sample of the substance you prepared or isolated. How this sample is labeled is very important. Again, learning a correct method of labeling bottles and vials can save time in the laboratory, because fewer mistakes will be made. More importantly, learning to label properly can decrease the danger inherent in having samples of material that cannot be identified correctly at a later date.

Solid materials should be stored and submitted in containers that permit the substance to be removed easily. For this reason, narrow-mouthed bottles or vials are not used for solid substances. Liquids should be stored in containers that will

not let them escape through leakage. Be careful not to store volatile liquids in containers that have plastic caps, unless the cap is lined with an inert material such as Teflon. Otherwise, the vapors from the liquid are likely to contact the plastic and dissolve some of it, thus contaminating the substance being stored.

On the label, print the name of the substance, its melting or boiling point, the actual and percentage yields, and your name. An illustration of a properly prepared label follows:

<p>Isopentyl Acetate BP 140°C Yield 3.81 g (42.2%) Joe Schmedlock</p>

Laboratory Glassware: Care and Cleaning

Because your glassware is expensive and you are responsible for it, you will want to give it proper care and respect. If you read this chapter carefully and follow the procedures presented here, you may be able to avoid some unnecessary expense. You may also save time because cleaning problems and replacing broken glassware are time-consuming.

If you are unfamiliar with the equipment found in an organic laboratory or are uncertain about how such equipment should be treated, this chapter provides some useful information. It includes topics such as cleaning glassware, caring for glassware when using corrosive or caustic reagents, and assembling components from your organic laboratory kit. At the end of this section are illustrations and names of most of the equipment you are likely to find in your drawer or locker.

3.1 Cleaning Glassware

Glassware can be cleaned easily if you clean it immediately. It is good practice to do your “dishwashing” right away. With time, organic tarry materials left in a container begin to attack the surface of the glass. The longer you wait to clean glassware, the more extensively this interaction will have progressed. Cleaning is then more difficult because water will no longer wet the surface of the glass as effectively. If you can’t wash your glassware immediately after use, soak the dirty pieces in soapy water. A half-gallon plastic container is convenient for soaking and washing glassware. Using a plastic container also helps prevent the loss of small pieces of equipment used in microscale techniques.

Various soaps and detergents are available for washing glassware. They should be tried first when washing dirty glassware. Organic solvents can also be used because the residue remaining in dirty glassware is likely to be soluble in some organic solvent. After the solvent has been used, the conical vial or flask probably will have to be washed with soap and water to remove the residual solvent. When you use solvents in cleaning glassware, use caution, because the solvents are hazardous (see Technique 1). Use fairly small amounts of a solvent for cleaning purposes. Usually 1–2 mL will be sufficient. Acetone is commonly used, but it is expensive. Your **wash acetone** can be used effectively several times before it is “spent.” Once your acetone is spent, dispose of it as your instructor directs. If acetone does not work, other organic solvents such as methylene chloride or toluene can be used.

CAUTION



Acetone is very flammable. Do not use it around flames.

For troublesome stains and residues that adhere to the glass despite your best efforts, use a mixture of sulfuric acid and nitric acid. Cautiously add about 20 drops of concentrated sulfuric acid and 5 drops of concentrated nitric acid to the flask or vial.

CAUTION

You must wear safety glasses when you are using a cleaning solution made from sulfuric acid and nitric acid. Do not allow the solution to come into contact with your skin or clothing. It will cause severe burns on your skin and create holes in your clothing. The acids may also react with the residue in the container.

Swirl the acid mixture in the container for a few minutes. If necessary, place the glassware in a warm-water bath and heat cautiously to accelerate the cleaning process. Continue heating until any sign of a reaction ceases. When the cleaning procedure is completed, decant the mixture into an appropriate waste container.

CAUTION

Do not pour the acid solution into a waste container that is intended for organic wastes.

Rinse the piece of glassware thoroughly with water and then wash with soap and water. For most common organic chemistry applications, any stains that survive this treatment are not likely to cause difficulty in subsequent laboratory procedures.

If the glassware is contaminated with stopcock grease (unlikely with the glassware recommended in this book), rinse the glassware with a small amount (1–2 mL) of methylene chloride. Discard the rinse solution into an appropriate waste container. Once the grease is removed, wash the glassware with soap or detergent and water.

3.2 Drying Glassware

The easiest way to dry glassware is to let it stand overnight. Store conical vials, flasks, and beakers upside down on a piece of paper towel to permit the water to drain from them. Drying ovens can be used to dry glassware if they are available, and if they are not being used for other purposes. Rapid drying can be achieved by rinsing the glassware with acetone and air-drying it or placing it in an oven. First, thoroughly drain the glassware of water. Then rinse it with one or two *small* portions (1–2 mL) of acetone. Do not use any more acetone than is suggested here. Return the used acetone to a waste acetone container for recycling. After you rinse the glassware with acetone, dry it by placing it in a drying oven for a few minutes or allow it to air-dry at room temperature. The acetone can also be removed by aspirator suction. In some laboratories, it may be possible to dry the glassware by blowing a *gentle* stream of dry air into the container. (Your laboratory instructor will indicate if you should do this.) Before drying the glassware with air, make sure that the air line is not filled with oil. Otherwise, the oil will be blown into the container, and you will have to clean it again. It is not necessary to blast the acetone out of the glassware with a wide-open stream of air; a gentle stream of air is just as effective and will not startle other people in the room.

Do not dry your glassware with a paper towel unless the towel is lint-free. Most paper will leave lint on the glass that can interfere with subsequent procedures performed in the equipment. Sometimes it is not necessary to dry a piece of equipment thoroughly. For example, if you are going to place water or an aqueous solution in a container, it does not need to be completely dry.

3.3 Ground-Glass Joints

It is likely that the glassware in your organic kit has **standard-taper ground-glass joints**. For example, the air condenser in the figure consists of an inner (male) ground-glass joint at the bottom and an outer (female) joint at the top. Each end

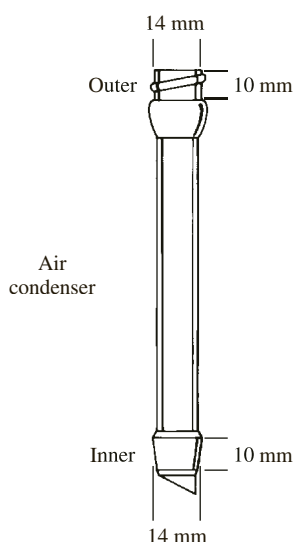


Figure 3.1
Illustration of F 14/10 inner and outer joints showing dimensions.

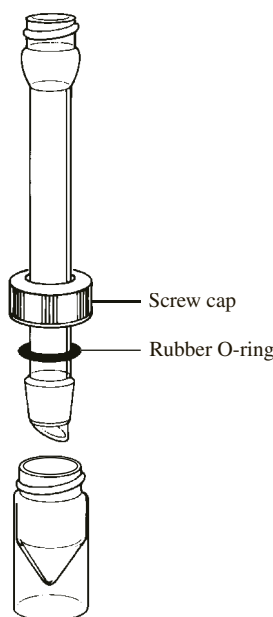


Figure 3.2
A microscale standard-taper joint assembly.

is ground to a precise size, which is designated by the symbol F followed by two numbers. A common joint size in microscale glassware is F 14/10. The first number indicates the diameter (in millimeters) of the joint at its widest point, and the second number refers to its length (see Figure 3.1). One advantage of standard-taper joints is that the pieces fit together snugly and form a good seal. In addition, standard-taper joints allow all glassware components with the same joint size to be connected, thus permitting the assembly of a wide variety of apparatus. One disadvantage of glassware with ground-glass joints, however, is that it is expensive.

Some pieces of glassware with ground-glass joints also have threads cast into the outside surface of the outer joints (see top of air condenser in Figure 3.1). The threaded joint allows the use of a plastic screw cap with a hole in the top to fasten two pieces of glassware together securely. The plastic cap is slipped over the inner joint of the upper piece of glassware, followed by a rubber O-ring (see Figure 3.2). The O-ring should be pushed down so that it fits snugly on top of the ground-glass joint. The inner ground-glass joint is then fitted into the outer joint of the bottom piece of glassware. The screw cap is tightened without excessive force to attach the entire apparatus firmly together. The O-ring provides an additional seal that makes this joint airtight. With this connecting system, it is unnecessary to use any type of grease to seal the joint. The O-ring *must be used* to obtain a good seal and to lessen the chances of breaking the glassware when you tighten the plastic cap.

It is important to make sure no solid or liquid is on the joint surfaces. Such material will lessen the efficiency of the seal, and the joints may leak. The presence of solid particles could cause the ground-glass joints to break when the plastic cap is tightened. Also, if the apparatus is to be heated, material caught between the joint surfaces will increase the tendency for the joints to stick. If the joint surfaces are coated with liquid or adhering solid, you should wipe them with a cloth or lint-free paper towel before assembling.

3.4 Separating Ground-Glass Joints

The most important thing you can do to prevent ground-glass joints from becoming “frozen,” or stuck together, is to disassemble the glassware as soon as possible after a procedure is completed. Even when this precaution is followed, ground-glass joints may become stuck tightly together. The same is true of glass stoppers in bottles or conical vials. Because microscale glassware is small and fragile, it is relatively easy to break a piece of glassware when trying to pull two pieces apart. If the pieces do not separate easily, you must be careful when you try to pull them apart. The best way is to hold the two pieces, with your hands touching, as close as possible to the joint. With a firm grasp, try to loosen the joint with a slight twisting motion (do not twist very hard). If this does not work, try to pull your hands apart without pushing sideways on the glassware.

If it is not possible to pull the pieces apart, the following methods may help. A frozen joint can sometimes be loosened if you tap it *gently* with the wooden handle of a spatula. Then, try to pull it apart as already described. If this procedure fails, you may try heating the joint in hot water or a steam bath. If this heating fails, the instructor may be able to advise you. As a last resort, you may try heating the joint in a flame. You should not try this unless the apparatus is hopelessly stuck, because heating by flame often causes the joint to expand rapidly and crack or break. If you use a flame, make sure the joint is clean and dry. Heat the outer part of the joint slowly, in the yellow portion of a low flame, until it expands and breaks away from the inner section. Heat the joint very slowly and carefully or it may break.

3.5 Etching Glassware Glassware that has been used for reactions involving strong bases such as sodium hydroxide or sodium alkoxides must be cleaned thoroughly *immediately* after use. If these caustic materials are allowed to remain in contact with the glass, they will etch the glass permanently. The etching makes later cleaning more difficult because dirt particles may become trapped within the microscopic surface irregularities of the etched glass. Furthermore, the glass is weakened, so the lifetime of the glassware is shortened. If caustic materials are allowed to come into contact with ground-glass joints without being removed promptly, the joints will become fused, or “frozen.” It is extremely difficult to separate fused joints without breaking them.

3.6 Assembling the Apparatus

Care must be taken when assembling the glass components into the desired apparatus. Always remember that Newtonian physics applies to chemical apparatus, and unsecured pieces of glassware are certain to respond to gravity. You should always clamp the glassware securely to a ring stand. Throughout this textbook, the illustrations of the various glassware arrangements include the clamps that attach the apparatus to a ring stand. You should assemble your apparatus using the clamps as shown in the illustrations.

3.7 Capping Conical Vials or Openings

The plastic screw caps used to join two pieces of glassware together can also be used to cap conical vials (see Figure 3.3) or other openings. A Teflon insert, or liner, fits inside the cap to cover the hole when the cap is used to seal a vial. Only one side of the liner is coated with Teflon. This side should always face toward the inside of the vial. (Note that the O-ring is not used when the cap is used to seal a vial.) To seal a vial, it is necessary to tighten the cap firmly but not too tightly. It is possible to crack the vial if you apply too much force. Some Teflon liners have a soft backing material (silicone rubber) that allows the liner to compress slightly when the cap is screwed down. It is easier to cap a vial securely with these liners without breaking the vial than with liners that have a harder backing material.

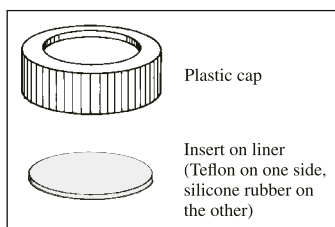



Figure 3.3
Plastic cap and Teflon insert.

3.8 Attaching Rubber Tubing to Equipment

When you attach rubber tubing to the glass apparatus or when you insert glass tubing into rubber stoppers, first lubricate the rubber tubing or the rubber stopper with either water or glycerin. Without such lubrication, it can be difficult to attach rubber tubing to the side arms of items of glassware such as condensers and filter flasks. Furthermore, glass tubing may break when it is inserted into rubber stoppers. Water is a good lubricant for most purposes. Do not use water as a lubricant when it might contaminate the reaction. Glycerin is a better lubricant than water and should be used when there is considerable friction between the glass and rubber. If glycerin is the lubricant, be careful not to use too much.


3.9 Description of Equipment

The components of the organic kit recommended for use in this textbook are given in Figures 3.4–3.7. Notice that most of the joints in these pieces of glassware are

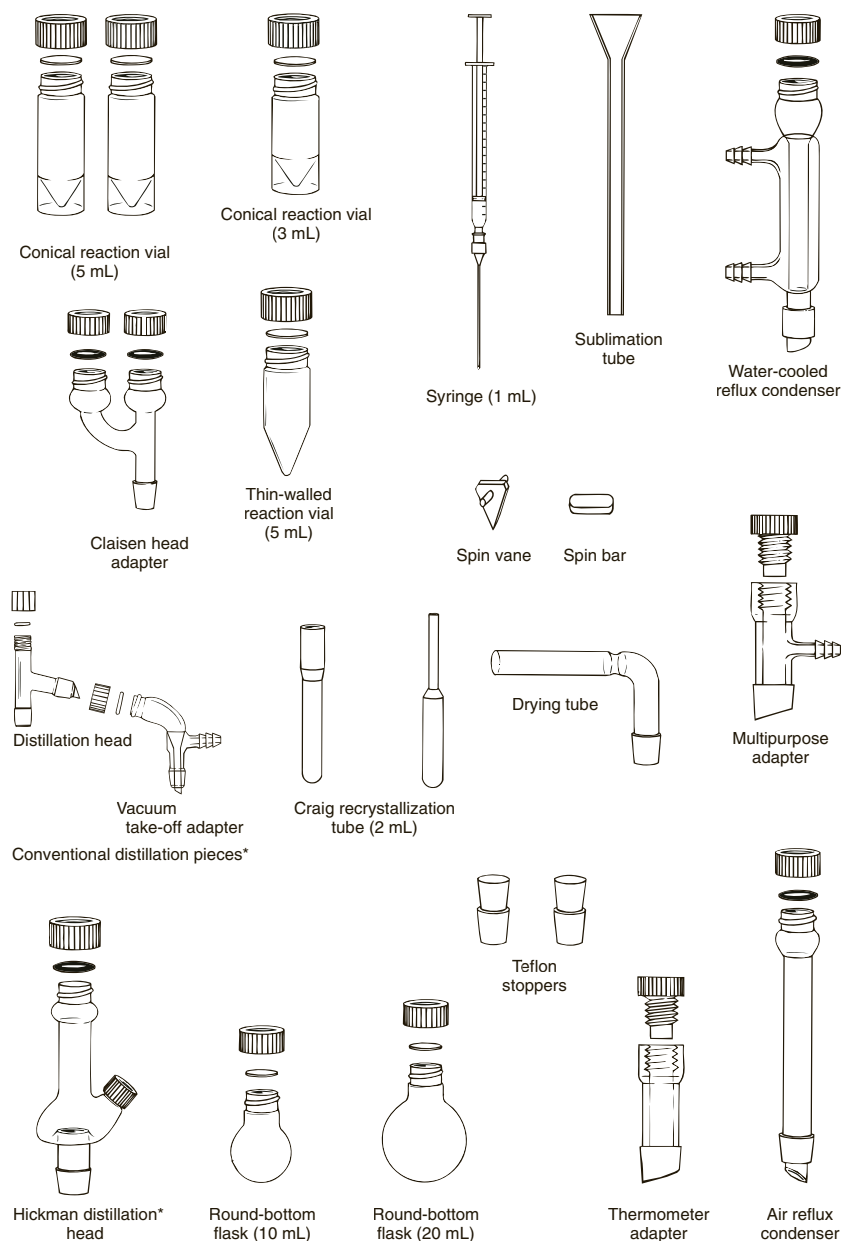

*Alternative types of distillation equipment are shown.

Figure 3.4
Components of a microscale organic kit.

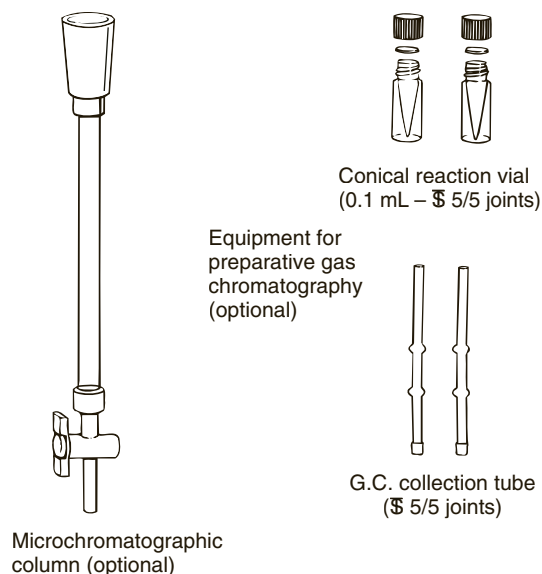**Figure 3.5**

Optional pieces of microscale glassware.

Note: The optional pieces of equipment shown in this figure are not part of the standard microscale kit. They must be purchased separately.

14/10, and all the outer joints are threaded. The organic kits used in your laboratory may have different joint sizes, or some of the outer joints may not be threaded. In particular, some older organic kits contain a number of pieces of glassware with 7/10 joints. These kits will work as well with the experiments in this book as the glassware recommended in the figures. In addition, there are microscale kits containing glassware that is connected without the use of ground-glass joints. The experiments in this book can also be performed with these glassware kits. Modifications with organic kits not containing the recommended glassware are discussed in the Technique chapters and in some of the experiments.

Figures 3.4–3.7 include glassware and equipment that are commonly used in the organic laboratory. Your glassware and equipment may vary slightly from the pieces shown on this spread and on the following pages.

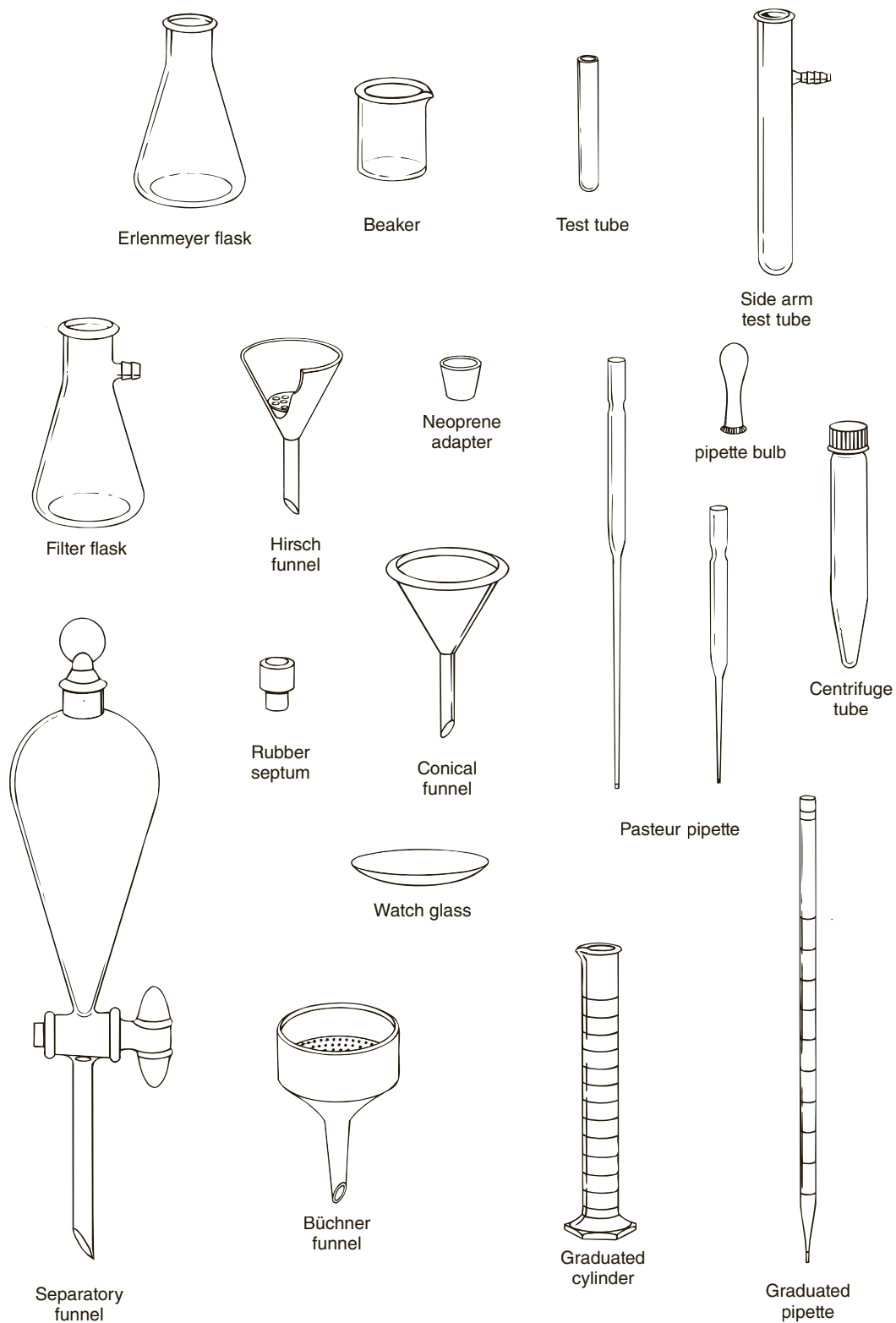


Figure 3.6
Equipment commonly used in the organic laboratory.

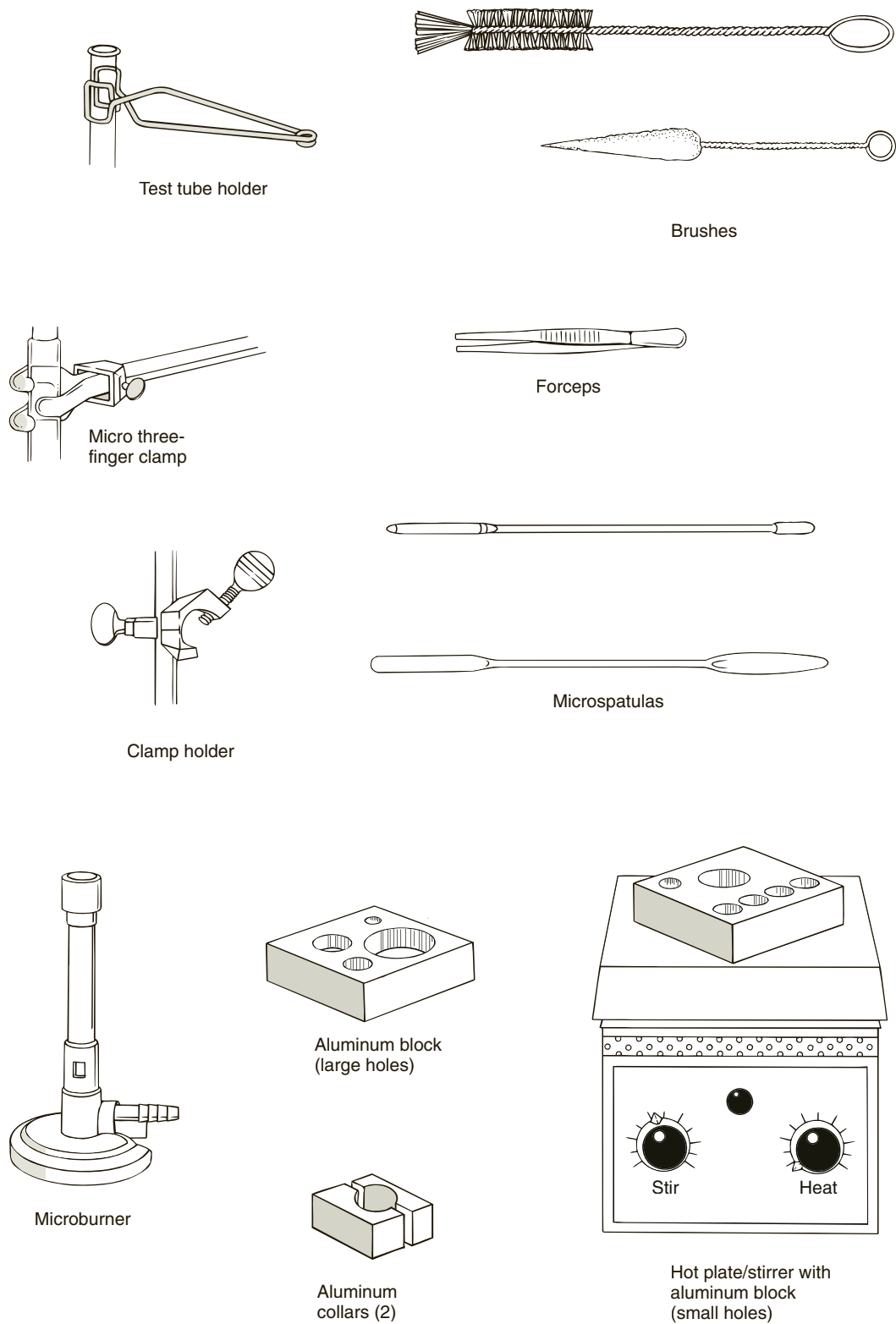


Figure 3.7
Equipment commonly used in the organic laboratory.

How to Find Data for Compounds: Handbooks and Catalogs

The best way to find information quickly on organic compounds is to consult a handbook. We will discuss the use of the *CRC Handbook of Chemistry and Physics*, *Lange's Handbook of Chemistry*, *The Merck Index*, and the *Aldrich Handbook of Fine Chemicals*. Complete citations to these handbooks are provided in Technique 29. Depending on the type of handbook consulted, the following information may be found:

- Name and common synonyms
- Formula
- Molecular weight
- Boiling point for a liquid or melting point for a solid
- Beilstein reference
- Solubility data
- Density
- Refractive index
- Flash point
- Chemical Abstracts Service (CAS) Registry Number
- Toxicity data
- Uses and synthesis

We often make use of the Internet to obtain information rapidly. To search the Internet, you can perform an online search using a browser (e.g., Mozilla Firefox, Internet Explorer, or Google Chrome) and a search engine (such as Google or Bing) to find the structures of compounds. Although this works well for obtaining the structures of organic compounds, some of the data, such as melting points and boiling points, may not be as reliable as data obtained from the handbooks listed above. If you are using the Internet to obtain data, make sure that you check several sources and confirm these values with those found in a handbook.

4.1 CRC Handbook of Chemistry and Physics

This is the handbook that is most often consulted for data on organic compounds. Although a new edition of the handbook is published each year, the changes that are made are often minor. An older copy of the handbook will often suffice for most purposes. In addition to the extensive tables of properties of organic compounds, the *CRC Handbook* includes sections on nomenclature and ring structures, an index of synonyms, and an index of molecular formulas.

Table 4.1 Examples of Names of Compounds in the *CRC Handbook*

Name of Organic Compound	Location in <i>CRC Handbook</i>
1-Chloropentane	Pentane, 1-chloro-
1,4-Dichlorobenzene	Benzene, 1,4-dichloro-
4-Chlorotoluene	Benzene, 1-chloro-4-methyl-
Ethanoic acid	Acetic acid
<i>tert</i> -Butyl acetate (ethanoate)	Acetic acid, 1,1-dimethylethyl ester
Ethyl propanoate	Propanoic acid, ethyl ester
Isopentyl alcohol	1-Butanol, 3-methyl-
Isopentyl acetate (banana oil)	1-Butanol, 3-methyl-, acetate
Salicylic acid	Benzoic acid, 2-hydroxy-
Acetylsalicylic acid (aspirin)	Benzoic acid, 2-acetyloxy-

The nomenclature used in this book most closely follows the Chemical Abstracts system of naming organic compounds. This system differs, but only slightly, from standard IUPAC nomenclature. Table 4.1 lists some examples of how some commonly encountered compounds are named in this handbook. The first thing you will notice is that this handbook is not like a dictionary. Instead, you must first identify the *parent* name of the compound of interest. The parent names are found in alphabetical order. Once the parent name is identified and found, then you look for the particular substituent or substituents that may be attached to this parent.

For most compounds, it is easy to find what you are looking for as long as you know the parent name. Alcohols are, as expected, named by IUPAC nomenclature. Notice in Table 4.1 that the branched-chain alcohol, isopentyl alcohol, is listed as 1-butanol, 3-methyl.

Esters, amides, and acid halides are usually named as derivatives of the parent carboxylic acid. Thus, in Table 4.1, you find ethyl propanoate listed under the parent carboxylic acid, propanoic acid. If you have trouble finding a particular ester under the parent carboxylic acid, try looking under the alcohol part of the name. For example, isopentyl acetate is not listed under acetic acid, as expected, but instead is found under the alcohol part of the name (see Table 4.1). Fortunately, this handbook has a Synonym Index that nicely locates isopentyl acetate for you in the main part of the handbook.

Once you locate the compound by its name, you will find the following useful information:

CRC number	This is an identification number for the compound. You can use this number to find the molecular structure located elsewhere in the handbook. This is especially useful when the compound has a complicated structure.
Name and synonym	The Chemical Abstracts name and possible synonyms.
Mol. form.	Molecular formula for the compound.
Mol. wt.	Molecular weight.
CAS RN	Chemical Abstracts Service Registry Number. This number is useful for locating additional information on the compound in the primary chemical literature (see Technique 29, Section 29.11).
mp/°C	Melting point of the compound in degrees Celsius.

bp/°C	Boiling point of the compound in degrees Celsius. A number without a superscript indicates that the recorded boiling point was obtained at 760 mmHg pressure (atmospheric pressure). A number with a superscript indicates that the boiling point was obtained at reduced pressure. For example, for an entry of 234, 122 ¹⁶ would indicate that the compound boils at 234 °C at 760 mmHg and 122 °C at 16 mmHg pressure.
Den/g cm ⁻³	Density of a liquid. A superscript indicates the temperature in degrees Celsius at which the density was obtained.
n _D	Refractive index determined at a wavelength of 589 nm, the yellow line in a sodium lamp (D line). A superscript indicates the temperature at which the refractive index was obtained (see Technique 24).
Solubility	<u>Solubility classification</u> <u>Solvent abbreviations</u> 1 = insoluble ace = acetone 2 = slightly soluble bz = benzene 3 = soluble chl = chloroform 4 = very soluble EtOH = ethanol 5 = miscible eth = ether 6 = decomposes hx = hexane
Beil. ref.	Beilstein reference. An entry of 4-02-00-00157 would indicate that the compound is found in the 4th supplement in Volume 2, with no subvolume, on page 157 (see Technique 29, Section 29.10 for details on the use of Beilstein).
Merck No.	<i>Merck Index</i> number in the 11th edition of the handbook. These numbers change each time a new edition of <i>The Merck Index</i> is issued.

Examples of sample handbook entries for isopentyl alcohol (1-butanol, 3-methyl) and isopentyl acetate (1-butanol, 3-methyl, acetate) are shown in Table 4.2.

4.2 Lange's Handbook of Chemistry

This handbook tends not to be as available as the *CRC Handbook*, but it has some interesting differences and advantages. *Lange's Handbook* has synonyms listed at the bottom of each page, along with structures of more complicated molecules. The most noticeable difference is in how compounds are named. For many compounds, the system lists names as they would appear in a dictionary. Table 4.3 lists examples of how some commonly encountered compounds are named in this handbook. Most often, you do not need to identify the *parent* name. Unfortunately, *Lange's Handbook* frequently uses common names that are becoming obsolete. For example, propionate is used rather than propanoate. Nevertheless, this handbook

Table 4.2 Properties of Isopentyl Alcohol and Isopentyl Acetate as Listed in the *CRC Handbook*

No.	Name Synonym	Mol. Form. Mol. Wt.	CAS RN mp/°C	Merck No. bp/°C	Beil. Ref. den/g cm ⁻³	Solubility n _D
3627	1-Butanol, 3-methyl Isopentyl alcohol	C ₅ H ₁₂ O 88.15	123-51-3 -117.2	5081 131.1	4-01-00-01677 0.8104 ²⁰	ace 4; eth 4; EtOH 4 1.4053 ²⁰
3631	1-Butanol, 3-methyl, acetate Isopentyl acetate	C ₇ H ₁₄ O ₂ 130.19	123-92-2 -78.5	4993 142.5	4-02-00-00157 0.876 ¹⁵	H ₂ O 2; EtOH 5; eth 5; ace 3 1.4000 ²⁰

Table 4.3 Examples of Names of Compounds in *Lange's Handbook*

Name of Organic Compound	Location in <i>Lange's Handbook</i>
1-Chloropentane	1-Chloropentane
1,4-Dichlorobenzene	1,4-Dichlorobenzene
4-Chlorotoluene	4-Chlorotoluene
Ethanoic acid	Acetic acid
<i>tert</i> -Butyl acetate (ethanoate)	<i>tert</i> -Butyl acetate
Ethyl propanoate	Ethyl propionate
Isopentyl alcohol	3-Methyl-1-butanol
Isopentyl acetate (banana oil)	Isopentyl acetate
Salicylic acid	2-Hydroxybenzoic acid
Acetylsalicylic acid (aspirin)	Acetylsalicylic acid

often names compounds as a practicing organic chemist would tend to name them. Notice how easy it is to find the entries for isopentyl acetate and acetylsalicylic acid (aspirin) in this handbook.

Once you locate the compound by its name, you will find the following useful information:

Lange's number	This is an identification number for the compound.
Name	See examples in Table 4.3.
Formula	Structures are drawn out. If they are complicated, then the structures are shown at the bottom of the page.
Formula weight	Molecular weight of the compound.
Beilstein reference	An entry of 2, 132 would indicate that the compound is found in Volume 2 of the main work on page 132. An entry of 3 ² , 188 would indicate that the compound is found in Volume 3 of the second supplement on page 188 (see Technique 29, Section 29.10 for details on the use of <i>Beilstein</i>).
Density	Density is usually expressed in units of g/mL or g/cm ³ . A superscript indicates the temperature at which the density was measured. If the density is also subscripted, usually 4°, it indicates that the density was measured at a certain temperature relative to water at its maximum density, 4°C. Most of the time you can simply ignore the subscripts and superscripts.
Refractive index	A superscript indicates the temperature at which the refractive index was determined (see Technique 24).
Melting point	Melting point of the compound in degrees Celsius. When a "d" or "dec" appears with the melting point, it indicates that the compound decomposes at the melting point. When decomposition occurs, you will often observe a change in color of the solid.
Boiling point	Boiling point of the compound in degrees Celsius. A number without a superscript indicates that the recorded boiling point was obtained at 760 mmHg pressure (atmospheric pressure). A number with a superscript indicates that the boiling point was obtained at reduced pressure. For example, an entry of 102 ^{11 mm} would indicate that the compound boils at 102°C at 11 mmHg pressure.

Flash point	This number is the temperature in degrees Celsius at which the compound will ignite when heated in air and a spark is introduced into the vapor. There are a number of different methods that are used to measure this value, so this number varies considerably. It gives a crude indication of flammability. You may need this information when heating a substance with a hot plate. Hot plates can be a serious source of trouble because of the sparking action that can occur with switches and thermostats used in hot plates.
Solubility in 100 parts solvent	Parts by weight of a compound that can be dissolved in 100 parts by weight of solvent at room temperature. In some cases, the values given are expressed as the weight in grams that can be dissolved in 100 mL of solvent. This handbook is not consistent in describing solubility. Sometimes gram amounts are provided, but in other cases the description will be more vague, using terms such as <i>soluble</i> , <i>insoluble</i> , or <i>slightly soluble</i> .

Solvent abbreviations

acet = acetone
 bz = benzene
 chl = chloroform
 aq = water
 alc = ethanol
 eth = ether
 HOAc = acetic acid

Solubility characteristics

i = insoluble
 s = soluble
 sls = slightly soluble
 vs = very soluble
 misc = miscible

Examples of sample handbook entries for isopentyl alcohol (3-methyl-1-butanol) and isopentyl acetate are shown in Table 4.4.

4.3 The Merck Index

The Merck Index is a very useful book because it has additional information not found in the other two handbooks. This handbook, however, tends to emphasize medically related compounds, such as drugs and biological compounds, although it also lists many other common organic compounds. It is not revised each year; new editions are published in five- or six-year cycles. It does not contain all of the compounds listed in *Lange's Handbook* or the *CRC Handbook*. However, for the compounds listed, it provides a wealth of useful information. The handbook will provide you with some or all of the following data for each entry.

- Merck number, which changes each time a new edition is issued
- Name, including synonyms and stereochemical designation
- Molecular formula and structure
- Molecular weight
- Percentages of each of the elements in the compound
- Uses
- Source and synthesis, including references to the primary literature
- Optical rotation for chiral molecules
- Density, boiling point, and melting point

Solubility characteristics, including crystalline form

Pharmacology information

Toxicity data

One of the problems with looking up a compound in this handbook is trying to decide the name under which the compound will be listed. For example, isopentyl alcohol can also be named as 3-methyl-1-butanol or isoamyl alcohol. In the 12th edition of the handbook, it is listed under the name isopentyl alcohol (#5212) on page 886. Finding isopentyl acetate is an even more challenging task. It is located in the handbook under the name isoamyl acetate (#5125) on page 876. Often, it is easier to look up the name in the name index or to find it in the formula index.

The handbook has some useful appendices that include the CAS registry numbers, a biological activity index, a formula index, and a name index that also includes synonyms. When looking up a compound in one of the indexes, you need to remember that the numbers provided are compound numbers, rather than page numbers. There is also a very useful section on organic name reactions that includes references to the primary literature.

4.4 Aldrich Handbook of Fine Chemicals

The *Aldrich Handbook* is actually a catalog of chemicals sold by the Aldrich Chemical Company. The company includes in its catalog a large body of useful data on each compound that it sells. Because the catalog is reissued each year at no cost to the user, you should be able to find an old copy when the new one is issued. As you are mainly interested in the data on a particular compound and not the price, an old volume is perfectly fine. Isopentyl alcohol is listed as 3-methyl-1-butanol, and isopentyl acetate is listed as isoamyl acetate in the *Aldrich Handbook*. The following includes some of the properties and information listed for individual compounds.

Aldrich catalog number

Name: Aldrich uses a mixture of common and IUPAC names. It takes a bit of time to master the names. Fortunately, the catalog does a good job of cross-referencing compounds and has a very good molecular formula index.

CAS Registry Number

Structure

Synonym

Formula weight

Boiling point/melting point

Index of refraction

Density

Beilstein reference

Merck reference

Infrared spectrum reference to the Aldrich Library of FT-IR spectra

NMR spectrum reference to the Aldrich Library of ^{13}C and ^1H FT-NMR spectra

Literature references to the primary literature on the uses of the compound

Toxicity

Safety data and precautions

Flash point

Prices of chemicals

Table 4.4 Properties of 3-Methyl-1-Butanol and Isopentyl Acetate as Listed in *Lange's Handbook*

No.	Name	Formula	Formula Weight	Beilstein Reference	Density	Refractive Index	Melting Point	Boiling Point	Flash Point	Solubility in 100 Parts Solvent
m155	3-methyl-1-butanol	$(\text{CH}_3)_2\text{CHCH}_2\text{CH}_2\text{OH}$	88.15	1,392	0.8129^{15}_4	1.4085^{15}	-117.2	132.0	45	2 aq; misc alc, bz, chl, eth, HOAc
i80	Isopentyl acetate	$\text{CH}_3\text{COOCH}_2\text{CH}_2\text{CH}(\text{CH}_3)_2$	130.19	2,132	0.876^{15}_4	1.4007^{20}	-78.5	142.0	80	0.25 aq; misc alc, eth

4.5 Strategy for Finding Information: Summary

Most students and professors find *The Merck Index* and *Lange's Handbook* easier and more intuitive to use than the *CRC Handbook*. You can go directly to a compound without rearranging the name according to the parent or base name followed by its substituents. Another great source of information is the *Aldrich Handbook*, which contains those compounds that are easily available from a commercial source. Many compounds are found in the *Aldrich Handbook* that you may never find in any of the other handbooks. The Sigma-Aldrich Web site (<http://www.sigmaaldrich.com/>) allows you to search by name, synonym, and catalog number.

PROBLEMS

- Using *The Merck Index*, find and draw structures for the following compounds:
 - atropine
 - quinine
 - saccharin
 - benzo[a]pyrene (benzpyrene)
 - itaconic acid
 - adrenosterone
 - chrysanthemic acid (chrysanthemumic acid)
 - cholesterol
 - vitamin C (ascorbic acid)
- Find the melting points for the following compounds in the *CRC Handbook*, *Lange's Handbook*, or the *Aldrich Handbook*:
 - biphenyl
 - 4-bromobenzoic acid
 - 3-nitrophenol
- Find the boiling point for each compound in the references listed in problem 2:
 - octanoic acid at reduced pressure
 - 4-chloroacetophenone at atmosphere and reduced pressure
 - 2-methyl-2-heptanol
- Find the index of refraction n_D and density for the liquids listed in problem 3.
- Using the *Aldrich Handbook*, report the specific rotations for the enantiomers of camphor.
- Read the section on carbon tetrachloride in *The Merck Index* and list some of the health hazards for this compound.

Measurement of Volume and Weight

Special care must be taken when working with small amounts of liquid or solids. In the typical microscale experiment, a student will use from 10 to 1000 mg of a liquid or solid. Specially designed microscale equipment will be used for these small-scale reactions. You may not be used to working with such small quantities, but after a while you will adjust to “thinking small.”

Liquids to be used for an experiment will usually be found in small containers in a hood. For experiments in this book, an automatic pipette, dispensing pump, or calibrated pipette will be used for measuring the volume of a liquid. It is critical that **limiting reactants** be weighed for accuracy purposes. *Do not calculate the weight using densities!* Measurement of a small volume of a liquid is subject to a large experimental error when converted to a weight using the density of a liquid. To determine the weight of a liquid when dealing with limiting reactants, preweigh the container before adding the liquid to the container and then reweigh the container after adding the liquid. This gives an **exact weight** and avoids the experimental error involved in using densities to calculate weights when working with smaller amounts of a liquid. For **nonlimiting** liquid reactants, you may calculate the weight of the liquid from the volume you have delivered using the density of the liquid and the following equation:

$$\text{Weight (g)} = \text{density (g/mL)} \times \text{volume (mL)}$$

Solids may be found near the balance. When an accurate measurement is required, solids must be weighed on a balance that reads to the nearest milligram (0.001 g) or tenth of a milligram (0.0001 g). To weigh a solid, place your conical vial or round-bottom flask in a small beaker and take these with you to the balance. Place a piece of paper that has been folded once on the balance pan. The folded paper will enable you to pour the solid into the conical vial or flask without spilling. Use the larger of your two spatulas to aid the transfer of the solid to the paper. Never weigh directly into a conical vial or flask and never pour, dump, or shake a material from a bottle. While still at the balance, carefully transfer the solid from the paper to your vial or flask. The vial or flask should be in a beaker while you are transferring the solid. The beaker traps any material that fails to make it into the container. It also supports the vial or flask so that it does not fall over. It is not necessary to obtain the exact amount specified in the experimental procedure, and trying to be exact requires too much time at the balance. For example, if you obtained 0.140 g of a solid, rather than the 0.136 g specified in a procedure, you could use it, but the actual amount weighed should be recorded in your notebook. Use the amount you weighed to calculate the theoretical yield, if this solid is the limiting agent.

Careless dispensing of liquids and solids is a hazard in any laboratory. When reagents are spilled, you may be subjected to an unnecessary health or fire hazard. In addition, you may waste expensive chemicals, destroy balance pans and clothing, and damage the environment. Always clean up any spills immediately.

5.1 Automatic Pipettes When available, an automatic pipette increases the speed of transfer of liquids from reagent bottles. These pipettes are expensive and must be shared by the entire laboratory. A number of types of units are available commercially. We describe the use of the continuously adjustable automatic pipette. This type of pipette can be adjusted for any volume within its defined range using a three- or four-digit readout. Several types of adjustable automatic pipettes are shown in Figure 5.1. The typical laboratory may have several units available: one 10–100 μL (0.01–0.10 mL) pipette for smaller volumes, and two 100–1000 μL (0.10–1.00 mL) pipettes for larger volumes. Disposable tips are available for each of these units and are color coded: yellow and blue for the small and large units, respectively. The automatic pipette is accurate with aqueous solutions, but it is not as accurate with organic liquids.

In most cases, the instructor will adjust the pipette so that it will deliver the desired volume. It will be placed in a convenient location near the reagent bottle, usually in a hood, and students will reuse the tip. Your instructor will give directions for the correct use of the automatic pipette. Students must practice using the automatic pipette. Remember that the automatic pipette is expensive and must be handled carefully. To protect the unit, you must always use a tip on the end of the pipette. Liquid must be drawn only into this plastic tip and never up into the unit itself. If this happens, you should notify your laboratory instructor immediately. Keep the pipette upright and immerse the tip just below the surface of the liquid. Automatic pipettes should never be used with corrosive liquids, such as sulfuric acid or hydrochloric acid.

5.2 Dispensing Pumps Dispensing pumps may be used in place of automatic pipettes when larger amounts (more than 0.1 mL) of liquids are being dispensed in the laboratory. The pumps

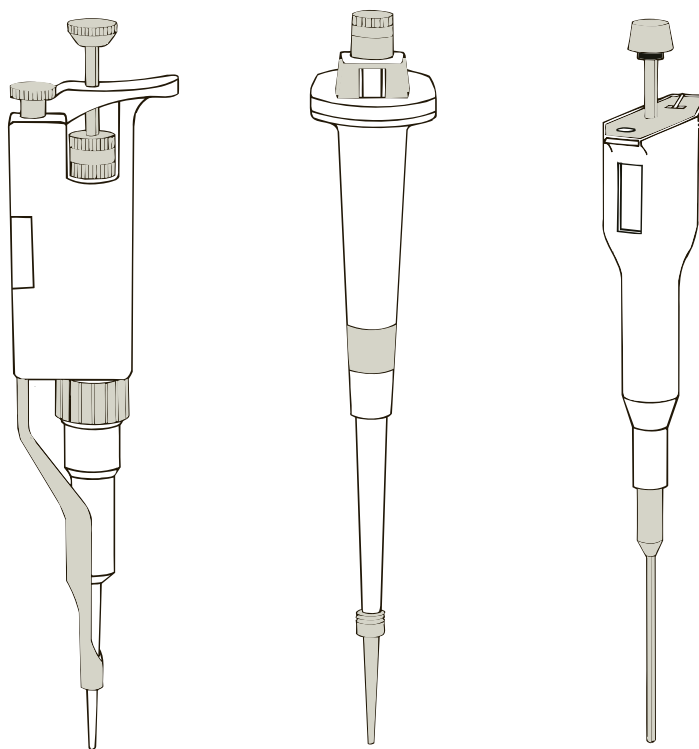


Figure 5.1
The adjustable automatic pipette.

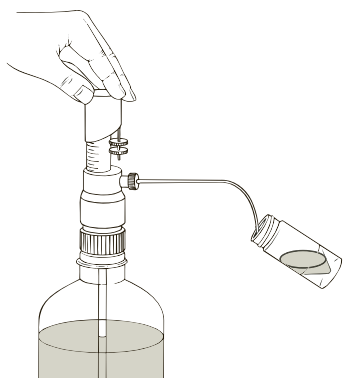


Figure 5.2
Dispensing pump.

are simple to operate, chemically inert, and accurate. Because the plunger assembly is made of Teflon, the dispensing pump may be used with most corrosive liquids and organic solvents. Dispensing pumps come in a variety of sizes, but the 1-, 2-, and 5-mL sizes are most useful in the microscale organic laboratory. The pump is attached to a bottle containing the liquid being dispensed. The liquid is drawn up from this reservoir into the pump assembly through a piece of inert plastic tubing.

Dispensing pumps are somewhat more difficult to adjust to the proper volume than automatic pipettes. Normally, the instructor or assistant will carefully adjust the unit to deliver the proper amount of liquid. As shown in Figure 5.2, the plunger is pulled up as far as it will travel to draw in the liquid from the glass reservoir. To expel the liquid from the spout into a container, slowly guide the plunger down. With low-viscosity liquids, the weight of the plunger will expel the liquid. With more viscous liquids, however, you may need to push the plunger gently to deliver the liquid into a container. Remove the last drop of liquid on the end of the spout by touching the tip on the interior wall of the container. When the liquid being transferred is a limiting reagent or when you need to know the weight precisely, you should weigh the liquid to determine the amount accurately.

As you pull up the plunger, look to see if the liquid is being drawn up into the pump unit. Some volatile liquids may not be drawn up in the expected manner, and you will observe an air bubble. Air bubbles are commonly observed when the pump has not been used for a while. The air bubble can be removed from the pump by dispensing and discarding several volumes of liquid to “reprime” the dispensing pump. Also check to see if the spout is filled completely with liquid. An accurate volume will not be dispensed unless the spout is filled with liquid before you lift up the plunger.

5.3 Graduated Pipettes

A suitable alternative to an automatic pipette or a dispensing pump is the graduated serological pipette. These *glass* pipettes are available commercially in a number of sizes. “Disposable” pipettes may be used many times and discarded only when the graduations become too faint to be seen. A good assortment of these pipettes consists of the following:

- 0.50-mL pipettes calibrated in 0.01-mL divisions (5/10 in 1/100 mL)
- 1.00-mL pipettes calibrated in 0.01-mL divisions (1 in 1/100 mL)
- 2.00-mL pipettes calibrated in 0.01-mL divisions (2 in 1/100 mL)

Liquids may be measured and transferred using a graduated pipette and a pipette pump. The style of pipette pump shown in Figure 5.3A is available in four sizes. The 2-mL size (blue) works well with the range of pipettes previously indicated. To fill the pipette, one simply rotates the knurled wheel forward so that the piston moves upward. The liquid is discharged by slowly turning the wheel backward until the proper amount of liquid has been expelled. The top of the pipette must be inserted securely into the pump and held there with one hand to obtain an adequate seal. The other hand is used to load and release the liquid.

The pipette pump shown in Figure 5.3B may also be used with graduated pipettes. The knob is turned counterclockwise to draw in the liquid, and then the liquid is released by turning the knob clockwise. With this style of pipette, the top of the pipette is held securely by a rubber O-ring, and it is easily handled with one hand. You should be certain that the pipette is held securely by the O-ring before using it. Disposable pipettes may not fit tightly in the O-ring, because they often have smaller diameters than nondisposable pipettes.

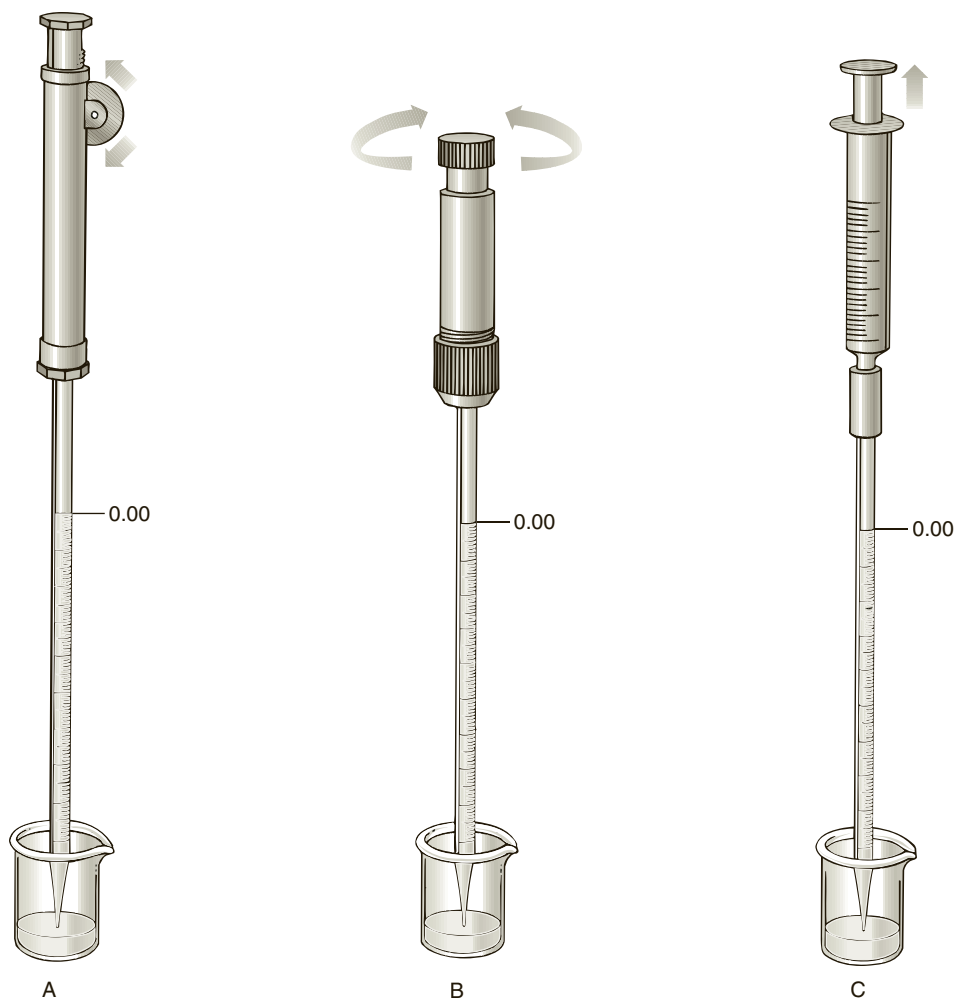


Figure 5.3
Pipette pumps.

A syringe may be used as a pipette pump, as shown in Figure 5.3C. In the design shown here, a 1- or 2-mL syringe is attached to the graduated pipette using a short piece of plastic tubing. The liquid is drawn up into the pipette when the plunger is pulled up, and it is expelled when the plunger is pushed down.

Excellent results may be obtained with graduated pipettes if you transfer by difference between marked calibrations and avoid transferring the entire contents of the pipette. When expelling the liquid, be sure to touch the tip of the pipette to the inside of the container before withdrawing the pipette. Graduated pipettes are commonly used when dispensing corrosive liquids, such as sulfuric acid or hydrochloric acid. The pipette will be supplied with a bulb or pipette pump.

Pipettes may be obtained in a number of styles, but only three types will be described here (Figure 5.4). One type of graduated pipette is calibrated “to deliver” (TD) its total capacity when the last drop is blown out. This style of pipette, shown in Figure 5.4A, is probably the most common type of graduated pipette in use in the laboratory; it is designated by two rings at the top. Of course, one does not need to transfer the entire volume to a container. To deliver a more accurate volume, you should transfer an amount less than the total capacity of the pipette using the graduations on the pipette as a guide.

Another type of graduated pipette is shown in Figure 5.4B. This pipette is calibrated to deliver its total capacity when the meniscus is located on the last

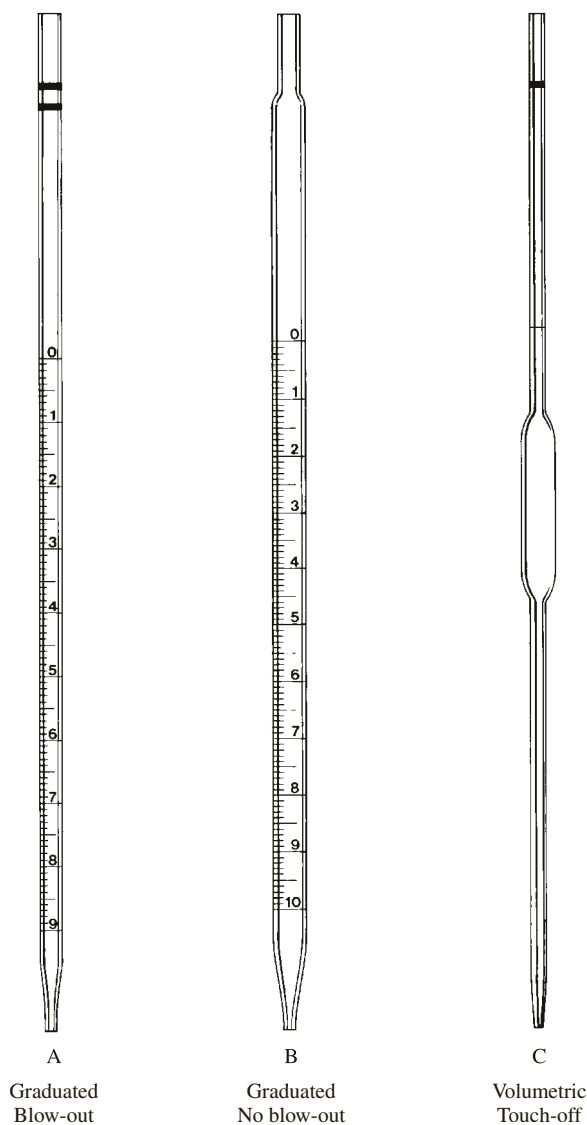


Figure 5.4
Pipettes.

graduation mark near the bottom of the pipette. For example, the pipette shown in Figure 5.4B delivers 10.0 mL of liquid when it has been drained to the point where the meniscus is located on the 10.0-mL mark. With this type of pipette, you must not drain the entire pipette or blow it out. In contrast, notice that the pipette shown in Figure 5.4A has its last graduation at 0.90 mL. The last 0.10-mL volume is blown out to give the 1.00-mL volume.

A nongraduated volumetric pipette is shown in Figure 5.4C. It is easily identified by the large bulb in the center of the pipette. This pipette is calibrated so that it will retain its last drop after the tip is touched on the side of the container. It must not be blown out. These pipettes often have a single colored band at the top that identifies it as a “touch-off” pipette. The color of the band is keyed to its total volume. This type of pipette is commonly used in analytical chemistry.

5.4 Pipettes The Pasteur pipette is shown in Figure 5.5A with a 2-mL rubber bulb attached. There are two sizes of pipettes: a long one (9 inch) and a short one (5¾ inch). It is important that the pipette bulb fit securely. You should not use a medicine dropper bulb, because of its small capacity. A Pasteur pipette is an indispensable piece of equipment for the routine transfer of liquids. It is also used for separations (Technique 12). Pasteur pipettes may be packed with cotton for use in gravity filtration (Technique 8) or packed with an adsorbent for small-scale column chromatography (Technique 19). Although they are considered disposable, you should be able to clean them for reuse as long as the tip remains unchipped.

A Pasteur pipette may be supplied by your instructor for dropwise addition of a particular reagent to a reaction mixture. For example, concentrated sulfuric acid is often dispensed in this way. When sulfuric acid is transferred,

take care to avoid getting the acid into the rubber or latex dropper bulb. It is best to avoid the rubber dropper bulb entirely by using one-piece transfer pipettes made entirely of polyethylene. These plastic pipettes are available in 1- or 2-mL sizes. They come from the manufacturers with approximate calibration marks stamped on them (Figure 5.5B).

Pipettes may be calibrated for use in operations where the volume does not need to be known precisely. Examples include measurement of solvents needed for extraction and for washing a solid obtained following crystallization. It is suggested that you calibrate several 5¾-inch pipettes following the procedure given in Experiment 1. A calibrated Pasteur pipette is shown in Figure 5.5C. Your instructor may provide you with a calibrated Pasteur pipette and bulb for transferring liquids where an accurate volume is not required. The pipette may be used to transfer a volume of 1.5 mL or less. You may find that the instructor has taped a test tube to the side of the storage bottle. The pipette is stored in the test tube with that particular reagent.

In general, Pasteur pipettes should not be used to measure volumes of reagents needed for organic reactions, because they are not accurate enough for this purpose.

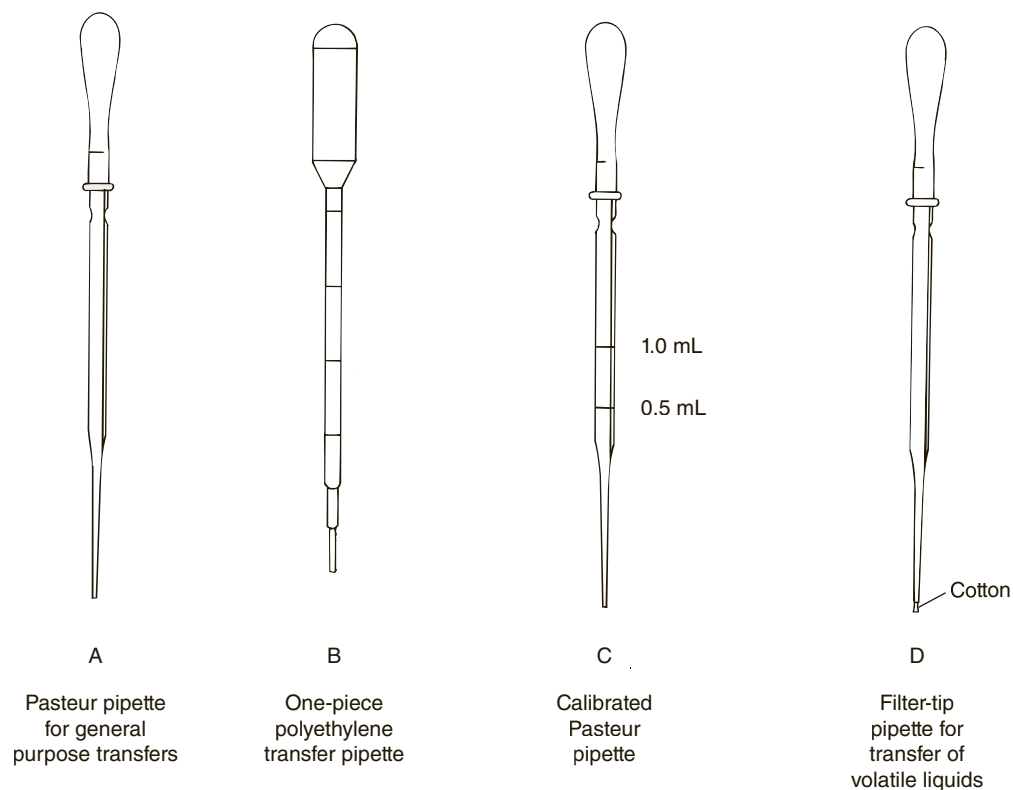


Figure 5.5
Pasteur and transfer pipettes.

In some cases, however, your instructor may have available a calibrated pipette for transferring nonlimiting reagents that may damage an automatic pipette. For example, a calibrated Pasteur pipette may be used with concentrated acids.

NOTE: You should not assume that a certain number of drops equals a 1-mL volume. The common rule that 20 drops equal 1 mL does not hold true for a Pasteur pipette!

A Pasteur pipette may be packed with cotton to create a filter-tip pipette as shown in Figure 5.5D. This pipette is prepared by the instructions given in Technique 8, Section 8.6. Pipettes of this type are useful in transferring volatile solvents during extractions and in filtering small amounts of solid impurities from solutions.

5.5 Syringes

Syringes may be used to add a pure liquid or a solution to a reaction mixture. They are especially useful when anhydrous conditions must be maintained. The needle is inserted through a septum, and the liquid is added to the reaction mixture. Although syringes come in a number of sizes, we will use a 1-mL unit in this textbook. Caution should be used with disposable syringes because they often use solvent-soluble rubber gaskets on the plungers. A syringe should be cleaned carefully after each use by drawing acetone or another volatile solvent into it and expelling the solvent with the plunger. Repeat this procedure several times to clean the syringe thoroughly. Draw air through the barrel with an aspirator to dry the syringe.

Syringes are usually supplied with volume graduations inscribed on the barrel. Large-volume syringes are not accurate enough to be used for measuring liquids in microscale experiments. A small microliter syringe, however, such as that used in gas chromatography, delivers a precise volume.

5.6 Graduated Cylinders

Graduated cylinders are used to measure relatively large volumes of liquids where accuracy is not required. For example, you could use a 10-mL graduated cylinder to obtain about 2 mL of a solvent for a crystallization procedure. You should use an automatic pipette, dispensing pump, or graduated pipette for accurate transfer of liquids in microscale work. Use a *clean and dry* Pasteur pipette to transfer the liquid from the storage container into the graduated cylinder. Do not attempt to pour the liquid directly into the cylinder from the storage bottle or you may spill the fluid. Some instructors may want you to pour some of the liquid into a beaker first and then use a Pasteur pipette to transfer the liquid to a graduated cylinder. Remember that you should not take more than you need. Excess material should never be returned to the storage bottle. Unless you can convince someone else to take it, it must be poured into the appropriate waste container. You should be frugal in estimating amounts needed.

5.7 Measuring Volumes with Conical Vials, Beakers, and Erlenmeyer Flasks

Conical vials, beakers, and Erlenmeyer flasks all have graduations inscribed on them. Beakers and flasks can be used to give only a crude approximation of the volume. They are much less precise than graduated cylinders for measuring volume. In some cases, a conical vial may be used to estimate volumes. For example, the graduations are sufficiently accurate for measuring a solvent needed to wash a solid obtained on a Hirsch funnel after a crystallization. You should use an automatic pipette, dispensing pump, or graduated transfer pipette for accurate measurement of liquids.

5.8 Balances

Solids and some liquids will need to be weighed on a balance that reads to at least the nearest milligram (0.001 g). A top-loading balance (see Figure 5.6) works well if the balance pan is covered with a plastic draft shield. The shield has a flap that opens to allow access to the balance pan. An analytical balance (see Figure 5.7) may

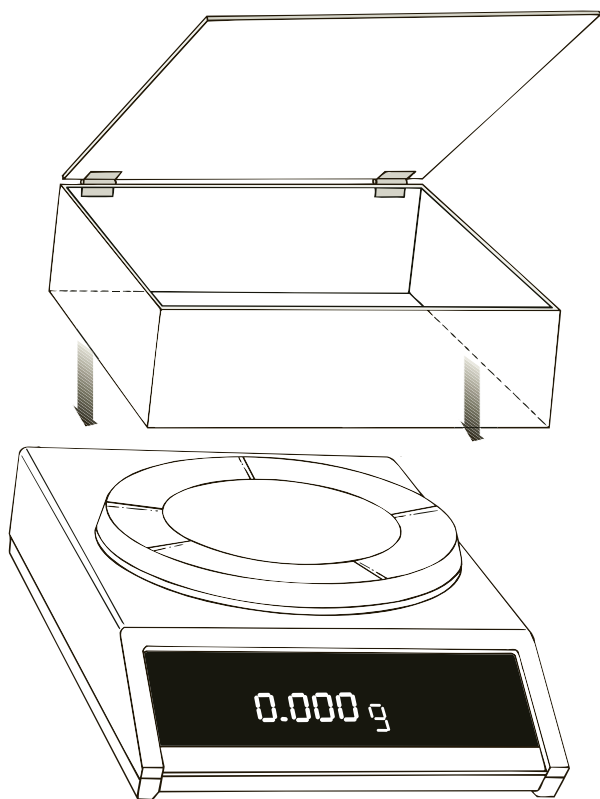


Figure 5.6
A top-loading balance with plastic draft shield.

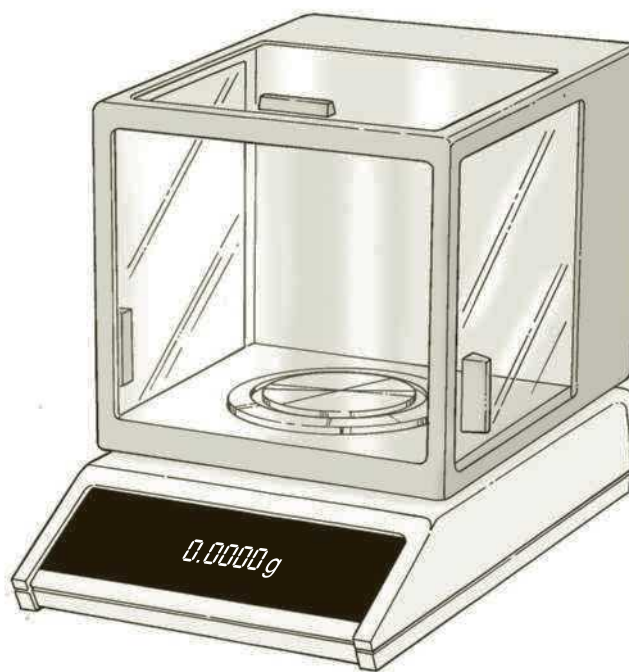


Figure 5.7
An analytical balance with glass draft shield.

also be used. This type of balance will weigh to the nearest tenth of a milligram (0.0001 g) when provided with a glass draft shield.

Modern electronic balances have a tare device that automatically subtracts the weight of a container or a piece of paper from the combined weight to give the weight of the sample. With solids, it is easy to place a piece of paper on the balance pan, press the tare device so that the paper appears to have zero weight, and then add your solid until the balance gives the weight you desire. You can then transfer the weighed solid to a container. You should always use a spatula to transfer a solid and never pour material from a bottle. In addition, solids must be weighed on paper and not directly on the balance pan. Remember to clean any spills.

With liquids, you should weigh the conical vial to determine the tare weight, transfer the liquid with an automatic pipette, dispensing pump, or graduated pipette into the vial, and then reweigh it. With liquids, it is usually necessary to weigh only the limiting reagent. The other liquids may be transferred using an automatic pipette, dispensing pump, or graduated pipette. Their weights can be calculated by knowing the volumes and densities of the liquids.


PROBLEMS

1. What measuring device would you use to measure the volume under each of the conditions described below? In some cases, there may be more than one answer to the question.
 - a. 5 mL of a solvent needed for a crystallization
 - b. 0.76 mL of a liquid needed for a reaction
 - c. 1 mL of a solvent needed for an extraction
2. Assume that the liquid used in part (b) is a limiting reagent for a reaction. What should you do after measuring the volume?
3. Calculate the weight of a 0.25-mL sample of each of the following liquids:
 - a. Diethyl ether (ether)
 - b. Methylene chloride (dichloromethane)
 - c. Acetone
4. A laboratory procedure calls for 0.146 g of acetic anhydride. Calculate the volume of this reagent needed in the reaction.
5. Criticize the following techniques:
 - a. A 100-mL graduated cylinder is used to measure accurately a volume of 2.8 mL.
 - b. A one-piece polyethylene transfer pipette (Figure 5.6B) is used to transfer precisely 0.75 mL of a liquid that is being used as the limiting reactant.
 - c. A calibrated Pasteur pipette (Figure 5.6C) is used to transfer 25 mL of a solvent.
 - d. The volume markings on a 100-mL beaker are used to transfer accurately 5 mL of a liquid.
 - e. An automatic pipette is used to transfer 10 mL of a liquid.
 - f. A graduated cylinder is used to transfer 0.126 mL of a liquid.
 - g. For a small-scale reaction, the weight of a liquid limiting reactant is calculated from its density and volume.

Heating and Cooling Methods


Most organic reaction mixtures need to be heated in order to complete the reaction. In general chemistry, you may have used a Bunsen burner for heating because nonflammable aqueous solutions were used. In an organic laboratory, however, the student must heat nonaqueous solutions that may contain *highly flammable* solvents. You *should not heat organic mixtures with a Bunsen burner* unless you are directed by your laboratory instructor. Open flames present a fire hazard. Whenever possible you should use one of the alternative heating methods, as described in the following sections.

6.1 Aluminum Block with Hot Plate/Stirrer

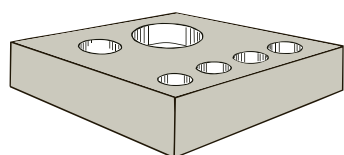

A

Aluminum block with small holes to fit Craig tube and 3-mL and 5-mL conical vials.

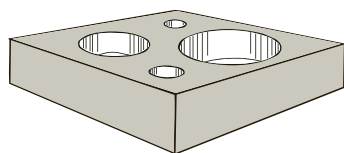

B

Aluminum block with large holes to fit 10-mL and 25-mL flasks

Most microscale organic laboratories now use an aluminum block and a hot plate, rather than a sand bath, for heating conical vials or flasks. There are several advantages to heating with an aluminum block. First, the metal will heat faster than a sand bath. Second, you can obtain a higher temperature with an aluminum block. Higher temperatures are often needed when distilling liquids with high boiling points at atmospheric pressure or under vacuum. Third, you can cool the aluminum rapidly by removing it with crucible tongs and immersing it in cold water.

Aluminum heating blocks can be fabricated readily in a machine shop or purchased from commercial suppliers.¹ The two aluminum blocks shown in Figure 6.1 will handle most heating applications in the laboratory. The block with the smaller holes will hold conical vials (Figure 6.1A). Holes have been drilled in the block so that different-sized conical vials will fit into the holes. This aluminum block may also be used in crystallizations using a Craig tube (Techniques 8 and 11). A hole is often provided for a mercury thermometer, but we do not recommend using it (see the caution box that follows). The aluminum block with the larger holes, as shown in Figure 6.1B, is designed to hold 10-, 20-, or 25-mL round-bottom flasks, as well as a thermometer.

Figure 6.2 shows a reaction mixture being heated with an aluminum block on a hot plate/stirrer unit. Also shown in Figure 6.2 is a split aluminum collar that may be used when very high temperatures are required. The collar is split to facilitate easy placement around a 5-mL conical vial. The collar helps distribute heat farther up the wall of the vial.

Because some hot plates vary widely in the temperature achieved for a given dial setting, some instructors may ask you to calibrate the hot plate so you have an approximate idea where to set the control on the hot plate to achieve a desired temperature. Place an aluminum block on the hot plate and

Figure 6.1
Aluminum heating blocks.

¹ The use of solid aluminum heating devices was developed by Siegfried Lodwig at Centralia College, Centralia, WA: Lodwig, S. N. "The Use of Solid Aluminum Heat Transfer Devices in Organic Chemistry Laboratory Instruction and Research." *J. Chem. Educ.* 66 (1989): 77.

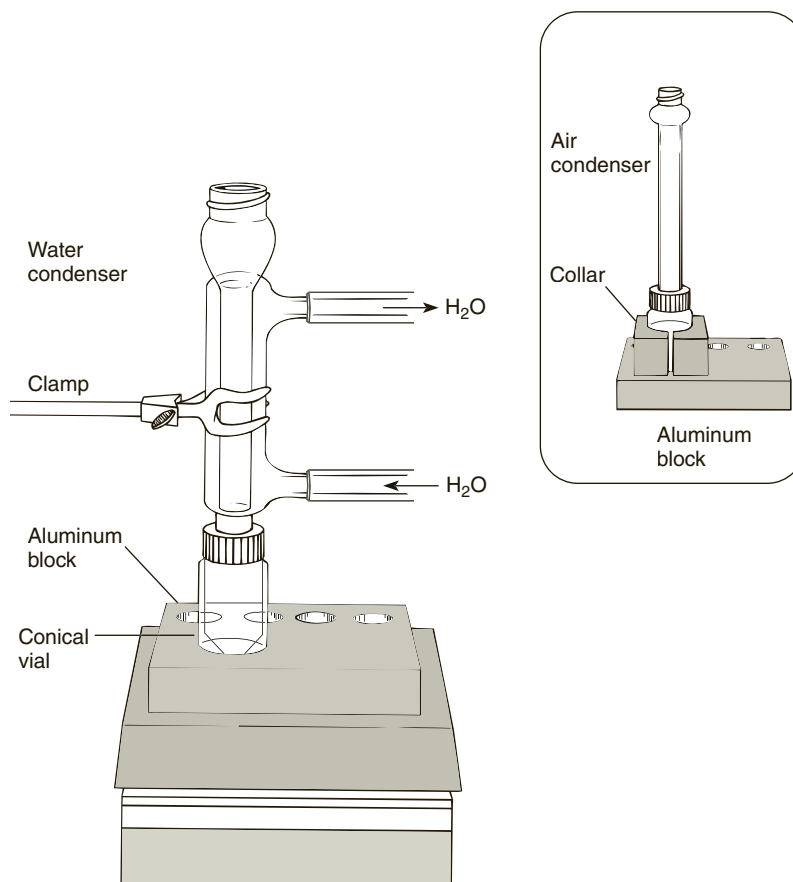


Figure 6.2
Heating with an aluminum block.

insert a non-mercury thermometer into the small hole in the block, as shown in Figure 6.2 (without the glassware). *Make sure the thermometer fits loosely in the hole or it may break.* Secure the thermometer with a clamp. Select five equally spaced settings on the heating control of the hot plate. Set the dial to the first of these settings and monitor the temperature recorded on the thermometer. When the thermometer reading arrives at a constant value, record this final temperature, along with the dial setting, in your notebook. Repeat this procedure with the remaining four settings and record the temperatures corresponding to the dial settings. Plot the data and keep it for future reference.

CAUTION



You should not use a mercury thermometer with an aluminum block. If it breaks, the mercury will vaporize on the hot surface. Instead, use a non-mercury glass thermometer, a metal dial thermometer, or a digital electronic temperature-measuring device.

To avoid the possibility of breaking a glass thermometer, your hot plate may have a hole drilled into the metal plate so that a metal dial thermometer can be inserted into the unit (Figure 6.3A). These metal thermometers, such as the one shown in Figure 6.3B, can be obtained in a number of temperature ranges. For example, a 0–250°C thermometer with 2-degree divisions can be obtained at a

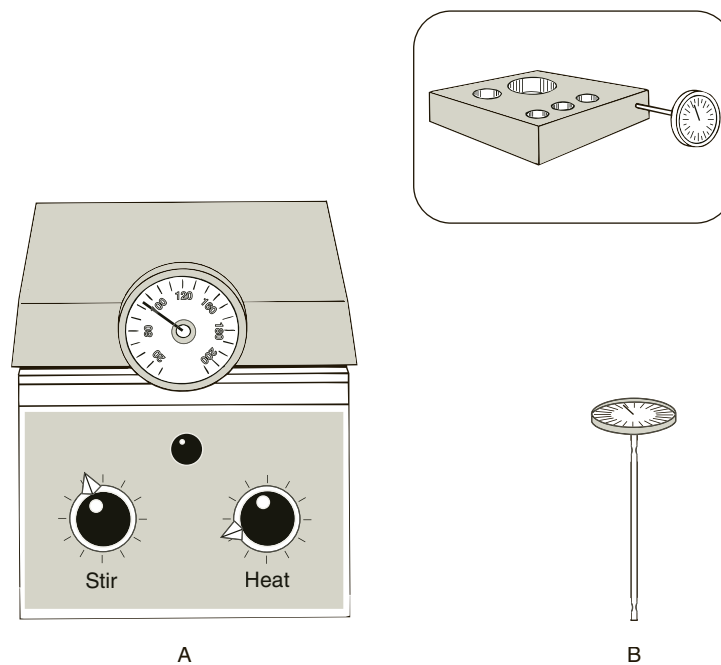


Figure 6.3
Dial thermometers.

reasonable price. Also shown in Figure 6.3 (inset) is an aluminum block with a small hole drilled into it so that a metal thermometer can be inserted.² An alternative to the metal thermometer is a digital electronic temperature-measuring device that can be inserted into the aluminum block or hot plate. It is strongly recommended that mercury thermometers be avoided when measuring the surface temperature of the hot plate or aluminum block. If a mercury thermometer is broken on a hot surface, you will introduce toxic mercury vapors into the laboratory. Non-mercury thermometers filled with high-boiling colored liquids are available as alternatives.

It is a good idea to use the same hot plate each time. It is likely that two hot plates of the same type may give different temperatures with an identical setting. Record the identification number printed on the unit that you are using in your notebook to ensure that you always use the same hot plate.

Although we may provide aluminum block temperatures for some experiments in this textbook, they should be taken as *approximate* values. You may need to adjust the temperature of the aluminum block appropriately to achieve the conditions you require. Each student must determine the actual temperature required to carry out a particular procedure. When a temperature is suggested, consider it as nothing more than a guide. Pay more attention to what is going on in your reaction vial or flask. If the temperature of your aluminum block equals the suggested temperature but the solution in the flask is not boiling (and you want it to boil), you clearly will need to increase the temperature of the aluminum block. Likewise, if the solution is boiling too rapidly, then you will need to reduce the temperature of the block.

When an aluminum block temperature is not given in the procedure and the liquid needs to be brought to a boil, you can determine the approximate

² Garner, C. M. "A Mercury-Free Alternative for Temperature Measurement in Aluminum Blocks." *J. Chem. Educ.* 68 (1991): A244.

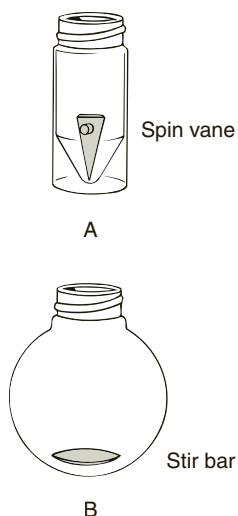


Figure 6.4
Methods of stirring
in a conical vial or
round-bottom flask.

6.2 Sand Bath with Hot Plate/Stirrer

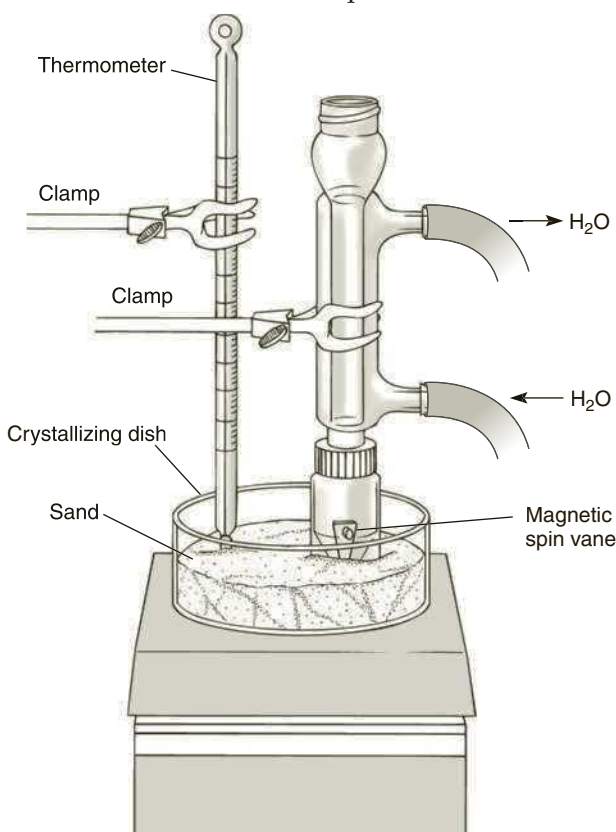


Figure 6.5
Heating with a sand bath.

setting from the boiling point of the liquid. Because the temperature inside the vial is lower than the aluminum block temperature, you should add at least 20°C to the boiling point of the liquid and set the aluminum block at this higher temperature. In fact, you may need to raise the temperature even higher than this value.

Many organic mixtures need to be stirred, as well as heated, to achieve satisfactory results. To stir a mixture, place a magnetic spin vane (Technique 7, Figure 7.8B) in a conical vial containing the reaction mixture as shown in Figure 6.4A. If the mixture is to be heated as well as stirred, attach a water condenser or an air condenser, as shown in Figure 6.2. With the combination stirrer/hot plate unit, it is possible to stir and heat a mixture simultaneously. Many reactions in this textbook are stirred continuously during the course of the reaction. With round-bottom flasks, a magnetic stir bar must be used to stir mixtures (Technique 7, Figure 7.8A). This is shown in Figure 6.4B. Many laboratories will have another aluminum block drilled to accommodate 10- and 25-mL round-bottom flasks. More uniform stirring will be obtained if the vial or flask is placed in the aluminum block so that it is centered on the hot plate. Mixing may also be achieved by boiling the mixture. A boiling stone (Technique 7, Section 7.4) should be added when a mixture is boiled without magnetic stirring.

The sand bath is used in some microscale laboratories to heat organic mixtures. Sand provides a clean way of distributing heat to a reaction mixture. To prepare a sand bath, place about a 1-cm depth of sand in a crystallizing dish or a Petri dish and then set the dish on a hot plate/stirrer unit. The apparatus is shown in Figure 6.5. Clamp the thermometer into position in the sand bath. You should calibrate the sand bath in a manner similar to that used with the aluminum block. Because sand heats more slowly than an aluminum block, you will need to begin heating the sand bath well before using it.

Do not heat the sand bath much above 200°C or you may break the dish. If you need to heat at very high temperatures, you should use an aluminum block rather than a sand bath (Section 6.1). With sand baths, it may be necessary to cover the dish with aluminum foil to achieve a temperature near 200°C. Keep in mind that the temperature obtained at a particular setting on the hot plate may vary for several reasons. First, you may place the thermometer at a different depth from time to time. Second, because of the relatively poor heat conduction of sand, you may obtain a different temperature in the conical vial, depending on the depth of the vial in the sand bath. Because of this poor heat conductivity, a temperature gradient is established within the sand bath. It is warmer near the bottom of the sand bath and cooler near the top for a given setting on the hot plate. To make use of this gradient, you may find it convenient to bury the vial or flask in the sand to heat a mixture more rapidly. Once the mixture is boiling, you can then slow the rate of heating by raising the vial or flask. These adjustments may be made easily and do not require a change in the setting on the hot plate.

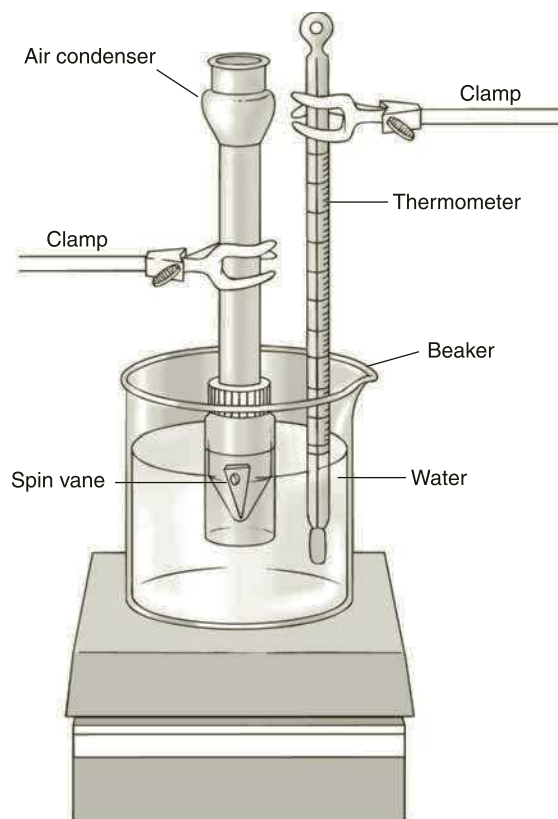


Figure 6.6
Water bath.

6.3 Water Bath with Hot Plate/Stirrer

The methods described previously may be used over a range of about 50°C to over 200°C. A hot water bath, however, may be a suitable alternative for temperatures below 80°C. A beaker (250 mL or 400 mL) is partially filled with water and heated on a hot plate. A thermometer is clamped into position in the water bath. You may need to cover the water bath with aluminum foil to prevent evaporation, especially at higher temperatures. The water bath is illustrated in Figure 6.6. A mixture can be stirred with a magnetic spin vane (Technique 7, Section 7.3). A hot water bath has some advantage over an aluminum block or a sand bath in that the temperature in the bath is uniform. In addition, it is sometimes easier to establish a lower temperature with a water bath than with other heating devices. Finally, the temperature of the reaction mixture will be closer to the temperature of the water, which allows for more precise control of the reaction conditions.

6.4 Flames

The simplest technique for heating mixtures is to use a Bunsen burner. Because of the high danger of fires, however, the use of the Bunsen burner should be strictly limited to those cases for which the danger of fire is low or for which no reasonable alternative source of heat is available. A flame should generally be used only to heat aqueous solutions or solutions with very high boiling points. You should always check with your instructor about using a burner. If you use a burner at your bench, great care should be taken to ensure that others in the vicinity are not using flammable solvents.

In heating a flask with a Bunsen burner, you will find that using a wire gauze can produce more even heating over a broader area. The wire gauze, when placed under the object being heated, spreads the flame to keep the flask from being heated in one small area only.

Bunsen burners may be used to prepare capillary micropipettes for thin-layer chromatography or to prepare other pieces of glassware requiring an open flame. For these purposes, burners should be used in designated areas in the laboratory and not at your laboratory bench.

6.5 Cold Baths

At times, you may need to cool a conical vial or flask below room temperature. A cold bath is used for this purpose. The most common cold bath is an **ice bath**, which is a highly convenient source of 0°C temperatures. An ice bath requires water along with ice to work well. If an ice bath is made up of only ice, it is not an efficient cooler because the large pieces of ice do not make good contact with the flask or vial. Enough water should be present with ice so that the flask is surrounded by water but not so much that the temperature is no longer maintained at 0°C . In addition, if too much water is present, the buoyancy of a flask resting in the ice bath may cause it to tip over. There should be enough ice in the bath to allow the flask to rest firmly.

For temperatures somewhat below 0°C , you may add some solid sodium chloride to the ice-water bath. The ionic salt lowers the freezing point of the ice so that temperatures from 0 to -10°C can be reached. The lowest temperatures are reached with ice-salt-water mixtures that contain relatively little water.

A temperature of -78.5°C can be obtained with solid carbon dioxide or dry ice. Large chunks of dry ice do not provide uniform contact with a flask being cooled. A liquid such as isopropyl alcohol is mixed with small pieces of dry ice to provide an efficient cooling mixture. Acetone and ethanol can be used in place of isopropyl alcohol. Be careful when handling dry ice because it can inflict severe frostbite. Extremely low temperatures can be obtained with liquid nitrogen (-195.8°C).

6.6 Steam Baths

The steam cone or steam bath is a good source of heat when temperatures less than 100°C are needed. Steam baths are used to heat reaction mixtures and solvents needed for crystallization. A steam cone and a portable steam bath are shown in Figure 6.7. These methods of heating have the disadvantage that water vapor may be introduced, through condensation of steam, into the mixture being heated. A slow flow of steam may minimize this difficulty.

Because water condenses in the steam line when it is not in use, it is necessary to purge the line of water before the steam will begin to flow. This purging should be accomplished before the flask is placed on the steam bath. The steam flow should be started with a high rate to purge the line; then the flow should be reduced to the desired rate. When using a portable steam bath, be certain that condensate (water)

is drained into a sink. Once the steam bath or cone is heated, a slow steam flow will maintain the temperature of the mixture being heated. There is no advantage to having a Vesuvius on your desk! An excessive steam flow may cause problems with condensation in the flask. This condensation problem can often be avoided by selecting the correct place at which to locate the flask on top of the steam bath.

The top of the steam bath consists of several flat concentric rings. The amount of heat

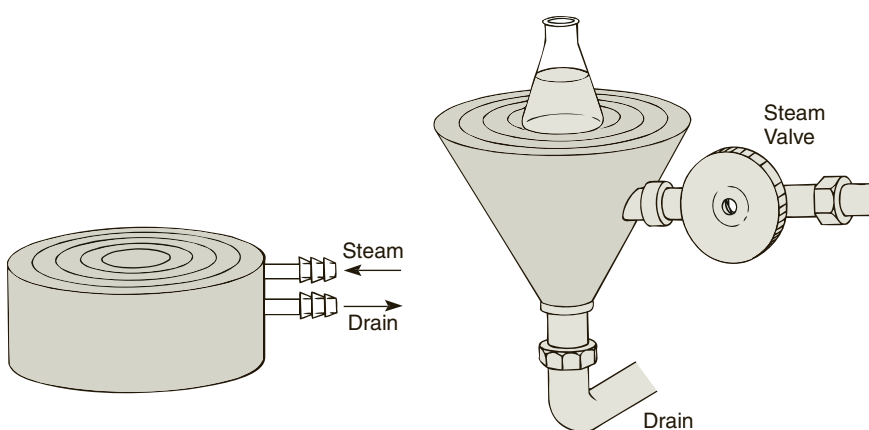


Figure 6.7
Steam bath and steam cone.

delivered to the flask being heated can be controlled by selecting the correct sizes of these rings. Heating is most efficient when the largest opening that will still support the flask is used. Heating large flasks on a steam bath while using the smallest opening leads to slow heating and wastes laboratory time.

PROBLEMS

1. What would be the preferred heating device(s) in each of the following situations?
 - a. Refluxing a solvent with a 56°C boiling point
 - b. Refluxing a solvent with a 120°C boiling point
 - c. Distilling a substance that boils at 220°C
2. Obtain the boiling points for the following compounds by using a handbook (Technique 4). In each case, suggest a heating device(s) that should be used for refluxing the substance.
 - a. Butyl benzoate
 - b. 1-Pentanol
 - c. 1-Chloropropane
3. What type of bath would you use to get a temperature of -10°C?
4. Obtain the melting point and boiling point for benzene and ammonia from a handbook (Technique 4) and answer the following questions.
 - a. A reaction was conducted in benzene as the solvent. Because the reaction was very exothermic, the mixture was cooled in an ice-salt-water bath. This was a bad choice. Why?
 - b. What bath should be used for a reaction that is conducted in *liquid* ammonia as the solvent?
5. Criticize the following techniques:
 - a. Refluxing a mixture that contains diethyl ether using a Bunsen burner
 - b. Refluxing a mixture that contains a large amount of toluene using a hot water bath
 - c. Using a mercury thermometer that is inserted into an aluminum block on a hot plate
 - d. Running a reaction with *tert*-butyl alcohol (2-methyl-2-propanol) that is cooled to 0°C in an ice bath

Reaction Methods

The successful completion of an organic reaction requires the chemist to be familiar with a variety of laboratory methods. These methods include operating safely, assembling the apparatus, heating and stirring reaction mixtures, adding liquid reagents, maintaining anhydrous and inert conditions in the reaction, and collecting gaseous products. Several techniques that are used in bringing a reaction to a successful conclusion are discussed here.

7.1 Assembling the Apparatus

Care must be taken when assembling the glass components into the desired apparatus. You should always remember that Newtonian physics applies to chemical apparatus, and unsecured pieces of glassware are certain to respond to gravity.

Assembling an apparatus in the correct manner requires that the individual pieces of glassware be connected to each other securely and the entire apparatus be held in the correct position. This can be accomplished by using **adjustable metal clamps** or a combination of adjustable metal clamps and **plastic joint clips**.

Two types of adjustable metal clamps are shown in Figure 7.1. Although these two types of clamps can usually be interchanged, the extension clamp is more commonly used to hold round-bottom flasks in place, and the three-finger clamp is frequently used to clamp condensers. Both types of clamps must be attached to a ring stand using a clamp holder, shown in Figure 7.1C.

A. Securing Macroscale Apparatus Assemblies

It is possible to assemble an apparatus using only adjustable metal clamps. An apparatus used to perform a distillation is shown in Figure 7.2. It is held together securely with three metal clamps. Because of the size of the apparatus and its geometry, the various clamps would likely be attached to three different ring stands. This apparatus would be somewhat difficult to assemble because it is necessary to ensure that the individual pieces stay together while securing and adjusting the clamps required to hold the entire apparatus in place. In addition, one must be careful not to bump any part of the apparatus or the ring stands after the apparatus is assembled.

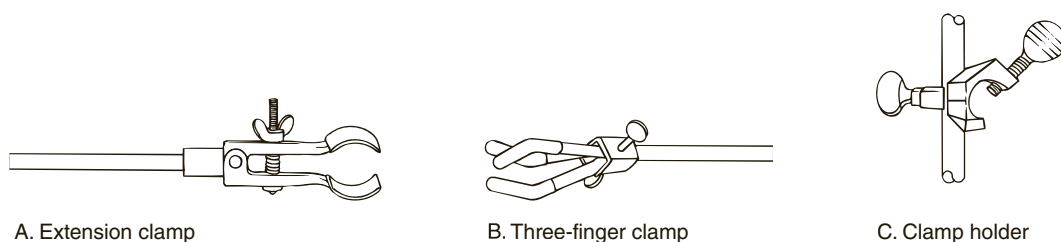


Figure 7.1
Adjustable metal clamps.


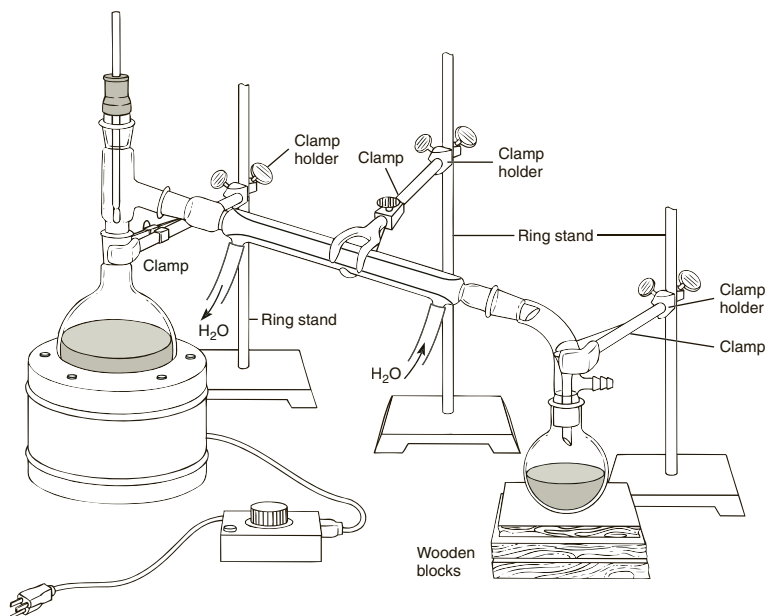


Figure 7.2
Distillation apparatus secured with metal clamps.

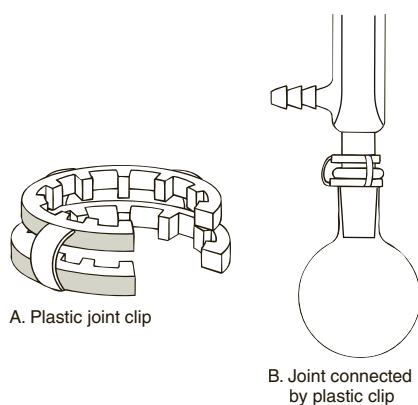

A. Plastic joint clip

B. Joint connected by plastic clip

Figure 7.3
Plastic joint clip.

A more convenient alternative is to use a combination of metal clamps and plastic joint clips. A plastic joint clip is shown in Figure 7.3A. These clips are easy to use (they just clip on), will withstand temperatures up to 140°C, and are durable. They hold together two pieces of glassware that are connected by ground-glass joints, as shown in Figure 7.3B. These clips come in different sizes to fit ground-glass joints of different sizes, and they are color coded for each size.

When used in combination with metal clamps, the plastic joint clips make it much easier to assemble most apparatus securely. There is less chance of dropping the glassware while assembling the apparatus, and once the apparatus is set up, it is more secure. Figure 7.4 shows the same distillation apparatus held in place with adjustable metal clamps and plastic joint clips.

To assemble this apparatus, first connect all of the individual pieces using the plastic clips. The entire apparatus is then connected to the ring stands using the adjustable metal clamps. Note that only two ring stands are required and the wooden blocks are not needed.

B. Securing Microscale Apparatus Assemblies

The glassware in most microscale kits is made with standard-taper ground joints. The most common joint size is $\text{F } 14/10$. Some microscale glassware with ground-glass joints also has threads cast into the outside surface of the outer joints (see top of air condenser in Figure 7.5). The threaded joint allows the use of a plastic screw cap with a hole in the top to fasten two pieces of glassware together securely. The plastic cap is slipped over the inner joint of the upper piece of glassware, followed by a rubber O-ring (see Figure 7.5). The O-ring should be pushed down so that it fits snugly on top of the ground-glass joint. The inner ground-glass joint is then fitted into the outer joint of the bottom piece of glassware. The screw cap is tightened, without excessive force, to attach the entire apparatus firmly together. The O-ring provides an additional seal that makes this joint airtight. With this connecting

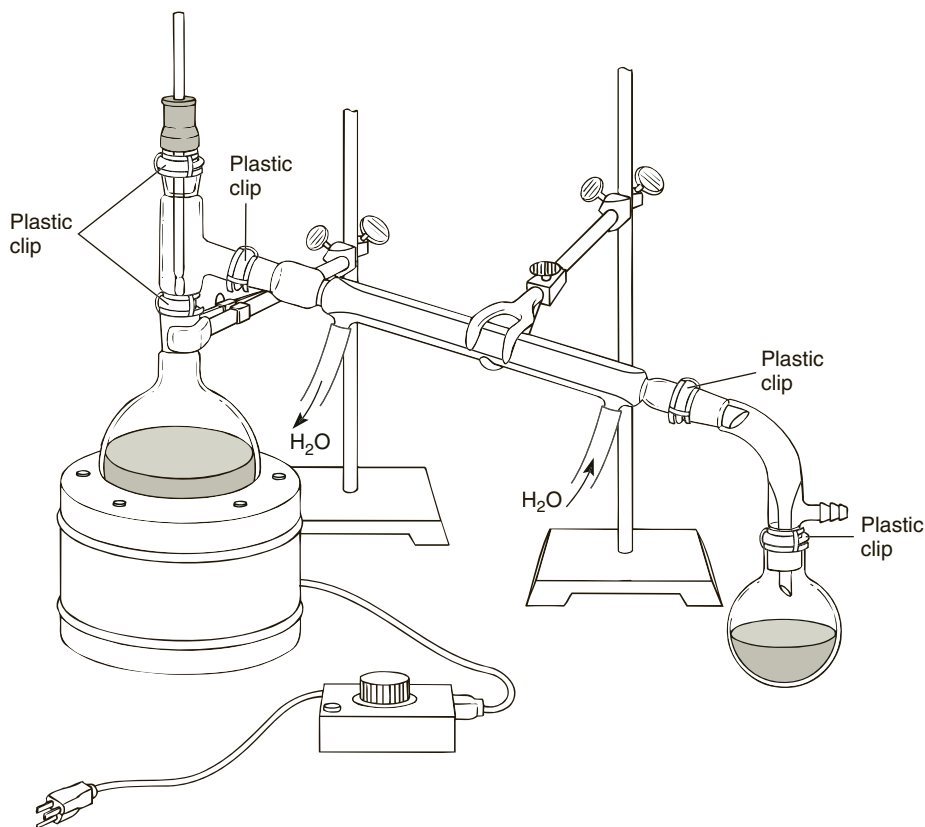


Figure 7.4
Distillation apparatus secured with metal clamps and plastic joint clips.

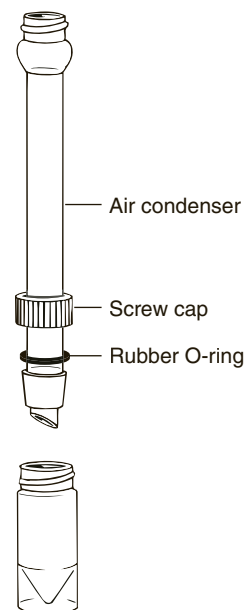


Figure 7.5
A microscale standard-taper joint assembly.

system, it is unnecessary to use any type of grease to seal the joint. The O-ring *must be used* to obtain a good seal and to lessen the chances of breaking the glassware when you tighten the plastic cap.

Microscale glassware connected in this fashion can be assembled easily. The entire apparatus is held together securely, and usually only one metal clamp is required to hold the apparatus onto a ring stand.

7.2 Heating under Reflux

Often we wish to heat a mixture for a long time and to leave it unattended. A **reflux apparatus** (see Figure 7.6) allows such heating. It also keeps the solvent from being lost by evaporation. A condenser is attached to the reaction vial or boiling flask.

Choice of Condenser. The condenser used in a reflux apparatus can be either of two types. An **air condenser** is simply a long tube. The surrounding air removes heat from the vapors within the tube and condenses them to liquid. A **water-jacketed condenser** consists of two concentric tubes with the outer cooling tube sealed onto the inner tube. The vapors rise within the inner tube, and water circulates through the outer tube. The circulating water removes heat from the vapors and condenses them. The air condenser is suitable for use with high-boiling liquids or with small quantities of material that are being heated gently. The water-jacketed condenser must be used when the vapors are difficult to condense, usually because the substance is volatile, or when vigorous boiling action is desired. In either case, the condenser prevents the vapors from escaping. Glassware assemblies using both air and water-jacketed condensers are shown in Figure 7.6A. The figure also shows a typical macroscale apparatus for heating large quantities of material under reflux (Figure 7.6B).

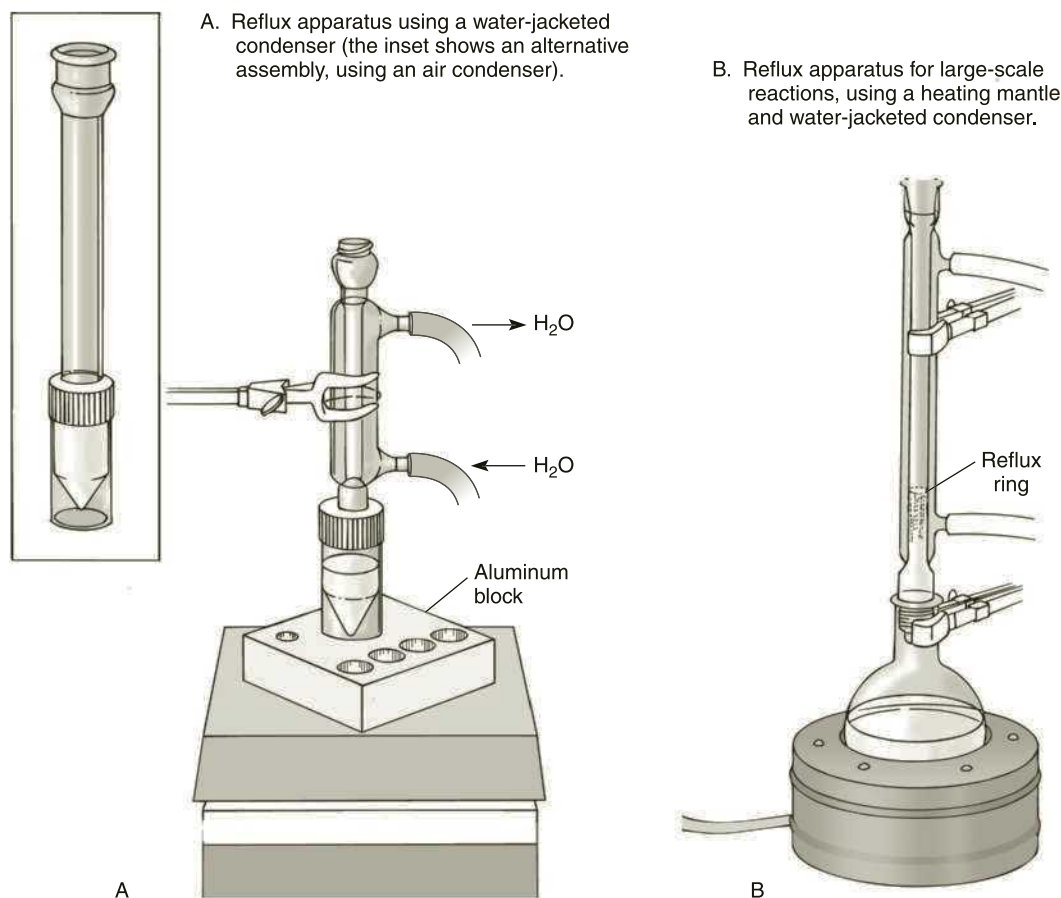


Figure 7.6
Heating under reflux.

When a water-jacketed condenser is used, the direction of the water flow should be such that the condenser will fill with cooling water. The water should enter the bottom of the condenser and leave from the top. The water should flow fast enough to withstand any changes in pressure in the water lines, but it should not flow any faster than absolutely necessary. An excessive flow rate greatly increases the chance of a flood, and high water pressure may force the hose from the condenser. Cooling water should be flowing before heating is begun! If the water is to remain flowing overnight, it is advisable to fasten the rubber tubing securely with wire to the condenser. If a flame is used as a source of heat, it is wise to use a wire gauze beneath the flask to provide an even distribution of heat from the flame. In most cases, an aluminum block, a sand bath, water bath, heating mantle, or steam bath is preferred over a flame.

Stirring. When heating a solution, always use a magnetic stirrer or a boiling stone (see Sections 7.3 and 7.4) to keep the solution from “bumping.”

Rate of Heating. If the heating rate has been correctly adjusted, the liquid being heated under reflux will travel only partway up the condenser tube before condensing. Below the condensation point, solvent will be seen running back into the flask; above it, the interior of the condenser will appear dry. The boundary between the two zones will be clearly demarcated, and a **reflux ring** or a ring of liquid will appear there. The reflux ring can be seen in Figure 7.6B. In heating under reflux,

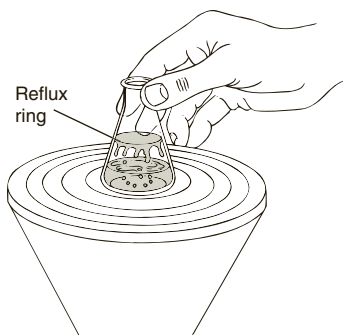


Figure 7.7
Tended reflux of small quantities on a steam cone (this can also be done with a hot plate).

the rate of heating should be adjusted so that the reflux ring is no higher than a third to half the distance to the top of the condenser. With microscale experiments, the quantities of vapor rising in the condenser frequently are so small that a clear reflux ring cannot be seen. In those cases, the heating rate must be adjusted so that the liquid boils smoothly but not so rapidly that solvent can escape the condenser. With such small volumes, the loss of even a small amount of solvent can affect the reaction. With large-scale reactions, the reflux ring is much easier to see, and one can adjust the heating rate more easily.

Tended Reflux. It is possible to heat small amounts of a solvent under reflux in an Erlenmeyer flask. With gentle heating, the evaporated solvent will condense in the relatively cold neck of the flask and return to the solution. This technique (see Figure 7.7) requires constant attention. The flask must be swirled frequently and removed from the heating source for a short period if the boiling becomes too vigorous. When heating is in progress, the reflux ring should not be allowed to rise into the neck of the flask.

How Do I Know How Hot to Heat It? A common problem that inexperienced students encounter when they assemble an apparatus for heating under reflux is that it is difficult to decide what temperature setting to use for heating the contents of a vial or flask to the desired temperature. This problem becomes more acute when the students attempt to reproduce the temperatures that are specified in the laboratory procedures of a textbook.

First, you should understand that the temperatures specified are only approximate suggestions. The actual temperature required to carry out a particular procedure must be determined for each student and each apparatus. When you see a temperature stipulated, consider it as nothing more than a guide. You will have to make adjustments to suit your own situation.

Second, you must always pay attention to what is going on in your reaction flask. If the temperature of your aluminum block or sand bath equals the suggested temperature, but the solution in your flask is not boiling, you clearly will have to increase the temperature of the heating device. Remember that what really matters is what is going on in the flask, not what the textbook says! The *external* temperature, as measured by a thermometer placed into the heating device, is not the important temperature. Far more critical is the temperature *inside* the flask, which may be considerably lower than the external temperature.

7.3 Stirring Methods

When a solution is heated, there is a danger that it may become superheated. When this happens, very large bubbles sometimes erupt violently from the solution; this is called **bumping**. Bumping must be avoided because of the risk that material may be lost from the apparatus, that a fire may start, or that the apparatus may break.

Magnetic stirrers are used to prevent bumping because they produce turbulence in the solution. The turbulence breaks up the large bubbles that form in boiling solutions. An additional purpose for using a magnetic stirrer is to stir the reaction to ensure that all the reagents are thoroughly mixed. A magnetic stirring system consists of a magnet that is rotated by an electric motor. The rate at which this magnet rotates can be adjusted by a potentiometric control. A small magnet, which is coated with a nonreactive material such as Teflon or glass, is placed in the flask. The magnet within the flask rotates in response to the rotating magnetic field caused by the motor-driven magnet. The result is that the inner magnet stirs the solution as it rotates. A common type of magnetic stirrer includes the stirring

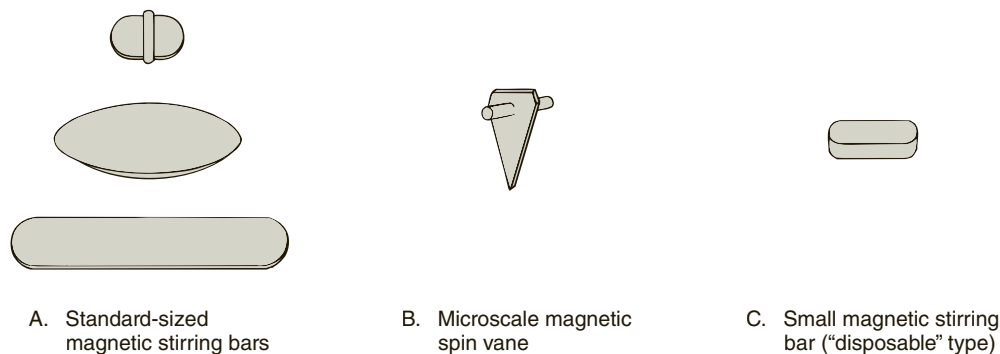


Figure 7.8
Magnetic stirring bars.

system within a hot plate. This type of hot plate/stirrer permits one to heat the reaction and stir it simultaneously. In order for the magnetic stirrer to be effective, the contents of the flask being stirred should be placed as close to the center of the hot plate as possible and not offset.

For macroscale apparatus, magnetic stirring bars of various sizes and shapes are available. For microscale apparatus, a **magnetic spin vane** is often used. It is designed to contain a tiny bar magnet and to have a shape that conforms to the conical bottom of a reaction vial. A small Teflon-coated magnetic stirring bar works well with very small round-bottom boiling flasks. Small stirring bars of this type (often sold as "disposable" stirring bars) can be obtained cheaply. A variety of magnetic stirring bars is illustrated in Figure 7.8.

There is also a variety of simple techniques that may be used to stir a liquid mixture in a centrifuge tube or conical vial. A thorough mixing of the components of a liquid can be achieved by repeatedly drawing the liquid into a Pasteur pipette and then ejecting the liquid back into the container by pressing sharply on the dropper bulb. Liquids can also be stirred effectively by placing the flattened end of a spatula into the container and twirling it rapidly.

7.4 Boiling Stones

A **boiling stone**, also known as a **boiling chip** or **Boileezer**, is a small lump of porous material that produces a steady stream of fine air bubbles when it is heated in a solvent. This stream of bubbles and the turbulence that accompanies it break up the large bubbles of gases in the liquid. In this way, it reduces the tendency of the liquid to become superheated, and it promotes the smooth boiling of the liquid. The boiling stone decreases the chances for bumping.

Two common types of boiling stones are carborundum and marble chips. Carborundum boiling stones are more inert, and the pieces are usually small, suitable for most applications. If available, carborundum boiling stones are preferred for most purposes. Marble chips may dissolve in strong acid solutions, and the pieces are larger. The advantage of marble chips is that they are cheaper.

Because boiling stones act to promote the smooth boiling of liquids, you should always make certain that a boiling stone has been placed in a liquid *before* heating is begun. If you wait until the liquid is hot, it may have become superheated. Adding a boiling stone to a superheated liquid will cause all the liquid to try to boil at once. The liquid, as a result, would erupt entirely out of the flask or froth violently.

As soon as boiling ceases in a liquid containing a boiling stone, the liquid is drawn into the pores of the boiling stone. When this happens, the boiling stone no longer can produce a fine stream of bubbles; it is spent. You may have to add a new boiling stone if you have allowed boiling to stop for a long period.

Wooden applicator sticks are used in some applications. They function in the same manner as boiling stones. Occasionally, glass beads are used. Their presence also causes sufficient turbulence in the liquid to prevent bumping.

7.5 Addition of Liquid Reagents

Liquid reagents and solutions are added to a reaction by several means, some of which are shown in Figure 7.9. For microscale experiments, the simplest approach is simply to add the liquid to the reaction by means of a Pasteur pipette. This method is shown in Figure 7.9A. In this technique, the system is open to the atmosphere. A second microscale method, shown in Figure 7.9B, is suitable for

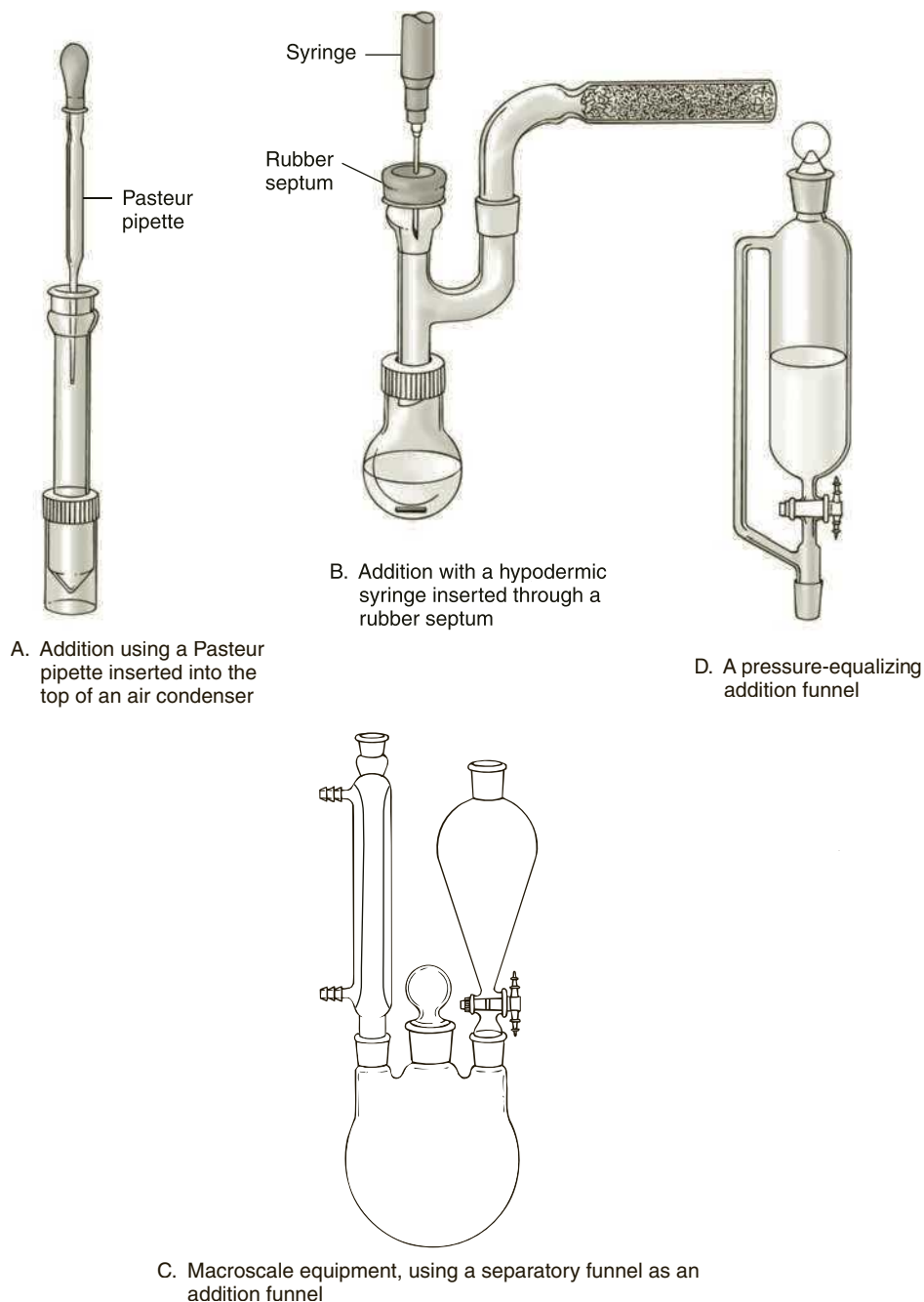


Figure 7.9
Methods of adding liquid reagents to a reaction.


experiments in which the reaction should be kept isolated from the atmosphere. In this approach, the liquid is kept in a hypodermic syringe. The syringe needle is inserted through a rubber septum, and the liquid is added dropwise from the syringe. The septum seals the apparatus from the atmosphere, which makes this technique useful for reactions that are conducted under an atmosphere of inert gas or where anhydrous conditions must be maintained. As an alternative, the rubber septum may be replaced by a cap and Teflon insert or liner. A disadvantage of the Teflon insert, however, is that the insert may no longer form an effective seal after being punctured by the needle.

The most common type of assembly for macroscale experiments is shown in Figure 7.9C. In this apparatus, a separatory funnel is attached to the side arm of a three-necked round-bottom flask. The separatory funnel must be equipped with a standard-taper, ground-glass joint to be used in this manner. The liquid is stored in the separatory funnel (which is called an **addition funnel** in this application) and is added to the reaction. The rate of addition is controlled by adjusting the stopcock. When it is being used as an addition funnel, the upper opening must be kept open to the atmosphere. If the upper hole is stoppered, a vacuum will develop in the funnel and will prevent the liquid from passing into the reaction vessel. Because the funnel is open to the atmosphere, there is a danger that atmospheric moisture can contaminate the liquid reagent as it is being added. To prevent this outcome, a drying tube (see Section 7.6) is attached to the upper opening of the addition funnel. The drying tube allows the funnel to maintain atmospheric pressure without allowing the passage of water vapor into the reaction.

Figure 7.9D shows an alternative type of addition funnel that is useful for reactions that must be maintained under an atmosphere of inert gas. This is the **pressure-equalizing addition funnel**. With this glassware, the upper opening is stoppered. The side arm allows the pressure above the liquid in the funnel to be in equilibrium with the pressure in the rest of the apparatus, and it lets the inert gas flow over the top of the liquid as it is being added.

7.6 Drying Tubes

With certain reactions, atmospheric moisture must be prevented from entering the reaction vessel. A **drying tube** can be used to maintain anhydrous conditions within the apparatus. Two types of drying tubes are shown in Figure 7.10. The typical drying tube is prepared by placing a small, loose plug of glass wool or cotton into the constriction at the end of the tube nearest the ground-glass joint

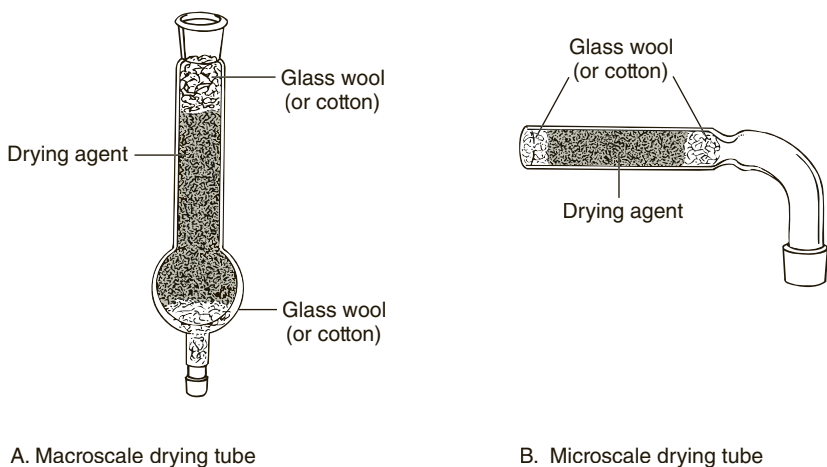

or hose connection. The plug is tamped gently with a glass rod or piece of wire to place it in the correct position. A drying agent, typically calcium sulfate (“Drierite”) or calcium chloride (see Technique 12, Section 12.9), is poured on top of the plug to the approximate depth shown in Figure 7.10. Another loose plug of glass wool or cotton is placed on top of the drying agent to prevent the solid material from falling out of the drying tube. The drying tube is then attached to the flask or condenser.

Air that enters the apparatus must pass through the drying tube.

Figure 7.10
Drying tubes.

7.7 Reactions Conducted under an Inert Atmosphere

The drying agent absorbs any moisture from air passing through it so that air entering the reaction vessel has had the water vapor removed from it.

Some reactions are sensitive to oxygen and water vapor present in air and require an inert atmosphere in order to obtain satisfactory results. The usual reactions in which it is desirable to exclude air often include organometallic reagents, such as organomagnesium or organolithium reagents, where water vapor and oxygen (air) react with these compounds. The most common inert gases available in a laboratory are nitrogen and argon, which are available in gas cylinders. Nitrogen is probably the gas most often used to carry out reactions under an inert atmosphere, although argon has a distinct advantage because it is denser than air. This allows the argon to push air away from the reaction mixture.

When laboratories are not equipped with individual gas lines to benches or hoods, it is useful to supply nitrogen or argon to the reaction apparatus using a balloon assembly (shown in Figure 7.11A). Your instructor will provide you with the apparatus.

Construct the balloon assembly by cutting off the top of a 3-mL disposable plastic syringe. Attach a small balloon snugly to the top of the syringe, securing it with a small rubber band that has been doubled to hold the balloon securely to the body of the syringe. Attach a needle to the syringe. Fill the balloon with the inert gas through the needle using a piece of rubber tubing attached to the gas source. When the balloon has been inflated to 2–3 inches in diameter, quickly pinch off the neck of the balloon while removing the gas source. Now push the needle into a rubber stopper to keep the balloon inflated. It is possible to keep an assembly like this filled with inert gas for several days without the balloon deflating.

Before you start the reaction, you may need to dry your apparatus thoroughly in an oven. Add all reagents carefully to avoid water. The following instructions are based on the assumption that you are using an apparatus consisting of a round-bottom flask equipped with a condenser. Attach a rubber septum to the top of your condenser. Now flush the air out of the apparatus with the inert gas. It is best not to use the balloon assembly for this purpose, unless you are using argon (see next paragraph). Instead, remove the round-bottom flask from the apparatus and, with the help of your instructor, flush it with the inert gas using a Pasteur pipette to bubble the gas through the solvent and reaction mixture in the flask. In this way, you can remove air from the reaction assembly before attaching the balloon assembly. Quickly reattach the flask to the apparatus. Pinch off the neck of the balloon between your fingers, remove the rubber stopper, and insert the needle into the rubber septum. The reaction apparatus is now ready for use.

When argon is employed as an inert gas, you can use the balloon assembly to remove air from the reaction apparatus in the following way. Insert the balloon assembly into the rubber septum as previously described. Also insert a second needle (no syringe attached) through the septum. The pressure from the balloon will force argon down the reflux condenser (argon is denser than air) and push the less dense air out through the second syringe needle. When the apparatus has been thoroughly flushed with argon, remove the second needle. Nitrogen does not work as well with this method because it is less dense than air, and it will be difficult to remove the air that is in contact with the reaction mixture in the round-bottom flask.

For reactions conducted at room temperature, you can remove the condenser shown in Figure 7.11A. Attach the rubber septum directly to the round-bottom flask and insert the needle of an argon-filled balloon

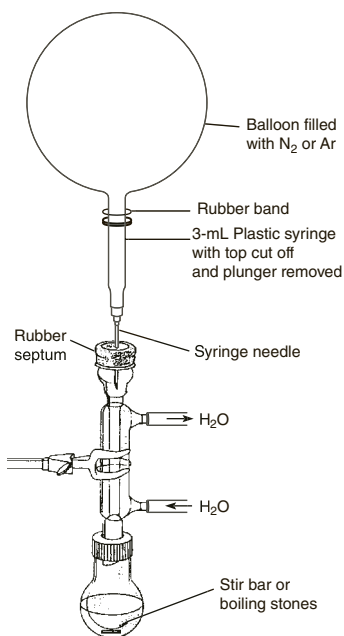


Figure 7.11A
Conducting a reaction under an inert atmosphere using a balloon assembly.

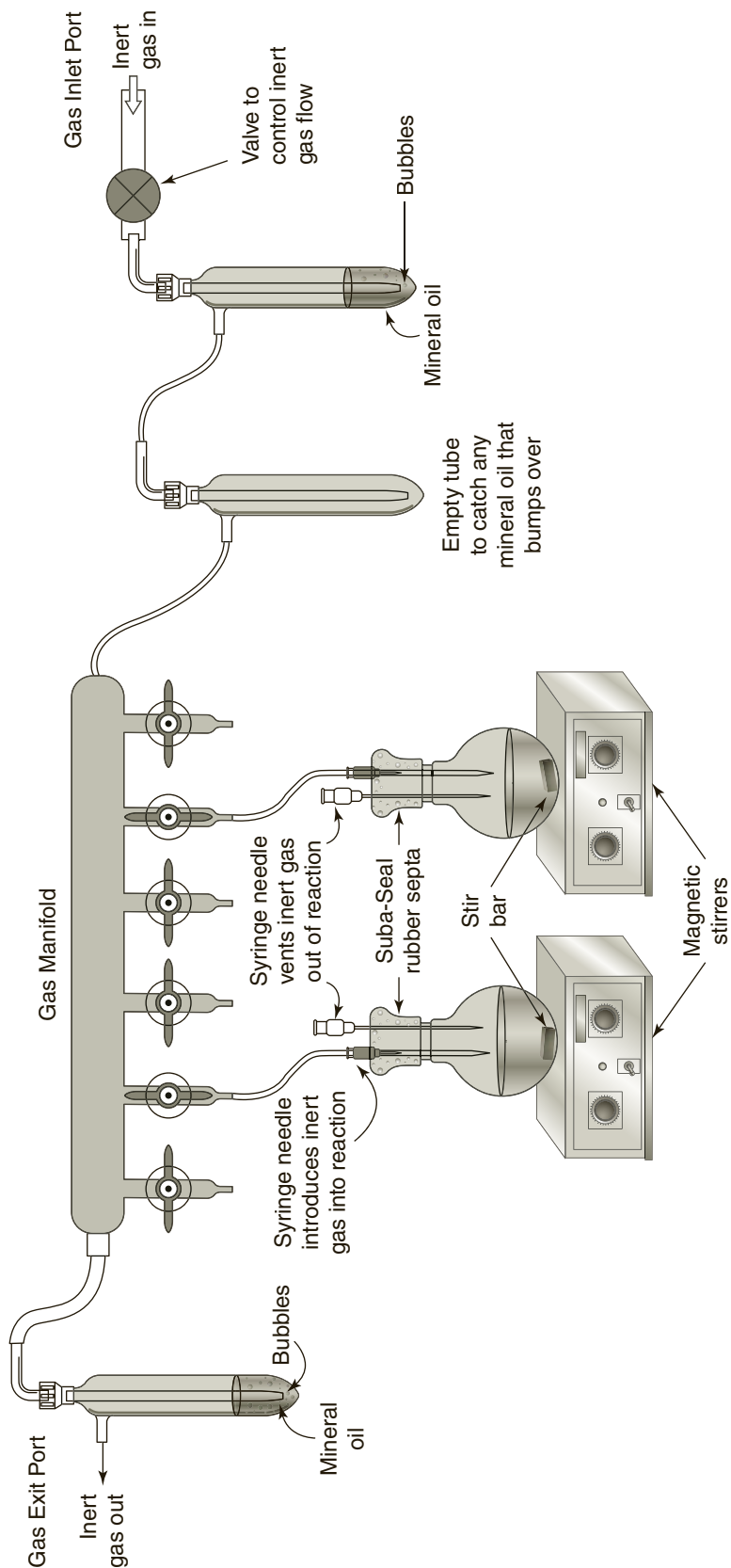


Figure 7.11B
Conducting a reaction under an inert atmosphere using a gas manifold.

assembly through the rubber septum. To flush the air out of the reaction flask, insert a second syringe needle into the rubber septum. Any air present in the flask will be flushed out through this second syringe needle, and the air will be replaced with argon. Now remove the second needle, and you have a reaction mixture free of air.

A gas manifold, shown in Figure 7.11B, provides another method of conducting reactions under an inert atmosphere. The gas manifold is equipped with multiple stopcocks so that more than one reaction can be conducted with a single tank of an inert gas, either nitrogen or argon. This setup allows a reaction to be maintained under anhydrous conditions in an oxygen-free environment. As shown in Figure 7.11B, the barrel of a cut-off syringe with an attached syringe needle is inserted into the rubber tubing. An inert gas is introduced through the Gas Inlet Port on the right side to expel air from the manifold. The mineral oil bubbler helps to monitor the gas flow into the manifold. The bubbling action in the inlet mineral oil bubbler must equal or exceed the bubbling action in the Gas Exit Port on the left side of Figure 7.11B.

Once the gas flow in the manifold is adjusted, one of the stopcocks on the manifold is opened, and the needle is inserted into a rubber septum on the reaction flask. The gas flow into the reaction flask is controlled with the stopcock on the manifold. Excess gas is vented through another syringe needle inserted into the rubber septum. Often the rate of bubbling observed in the Gas Exit Port will slow when the stopcock to the reaction flask is opened. The gas flow from the Gas Inlet Port will need to be increased in order to maintain bubbling in the Gas Exit Port. A complete description of the use of a gas manifold assembly is given in Experiment 63, *Synthesis of Naproxen by Palladium Catalysis*.

7.8 Capturing Noxious Gases

Many organic reactions involve the production of a noxious gaseous product. The gas may be corrosive, such as hydrogen chloride, hydrogen bromide, or sulfur dioxide, or it may be toxic, such as carbon monoxide. The safest way to avoid exposure to these gases is to conduct the reaction in a ventilated hood where the gases can be safely drawn away by the ventilation system.

In many instances, however, it is safe and efficient to conduct the experiment on the laboratory bench, away from the hood. This is particularly true when the gases are soluble in water. Some techniques for capturing noxious gases are presented in this section.

A. Drying Tube Method

Microscale experiments have the advantage that the amounts of gases produced are small. Hence, it is easy to trap them and prevent them from escaping into the laboratory. You can take advantage of the water solubility of corrosive gases such as hydrogen chloride, hydrogen bromide, and sulfur dioxide. A simple technique is to attach the drying tube (see Figure 7.10B) to the top of the reaction vial or condenser. The drying tube is filled with moistened glass wool. The moisture in the glass wool absorbs the gas, preventing its escape. To prepare this type of gas trap, fill the drying tube with glass wool and then add water dropwise to the glass wool until it has been moistened to the desired degree. Moistened cotton can also be used, although cotton will absorb so much water that it is easy to plug the drying tube.

When using glass wool in a drying tube, moisture from the glass wool must not be allowed to drain from the drying tube into the reaction. It is best to use a drying tube that has a constriction between the part where the glass wool is placed and the neck, where the joint is attached. The constriction acts as a partial barrier preventing the water from leaking into the neck of the drying tube. Make certain not to make the glass wool too moist.

B. External Gas Traps

Another approach to capturing gases is to prepare a trap that is separate from the reaction apparatus. The gases are carried from the reaction to the trap by means of tubing. There are several variations on this type of trap. One method that works well for microscale experiments is to place a thermometer adapter (Technique 14, Figure 14.9A) into the opening in the reaction apparatus. A Pasteur pipette is inserted upside down through the adapter, and a piece of fine flexible tubing is fitted over the narrow tip. It may be helpful to break the Pasteur pipette before using it for this purpose so that only the narrow tip and a short section of the barrel are used. The other end of the flexible tubing is placed through a large plug of moistened glass wool in a test tube. The water in the glass wool absorbs the water-soluble gases. This method is shown in Figure 7.12.

A variation on the Pasteur pipette method uses a hypodermic syringe needle inserted upside down (from the inside) through a rubber septum, which has been fitted over the opening at the top of the reaction apparatus. Flexible tubing, fitted over the syringe needle, leads to a trap such as the one using wet glass wool described previously. This variation is also shown in Figure 7.12.

Another alternative to the apparatus shown in Figure 7.12 is to use a multipurpose adapter in place of the thermometer adapter. The flexible tubing can be attached directly to the side arm of the multipurpose adapter, thus connecting the apparatus to the gas trap. If the multipurpose adapter is used for this purpose, the upper opening of the adapter must be closed; this is accomplished most easily by inserting a piece of glass rod or a short piece of glass tubing sealed at one end into the opening and tightening the fittings around it.

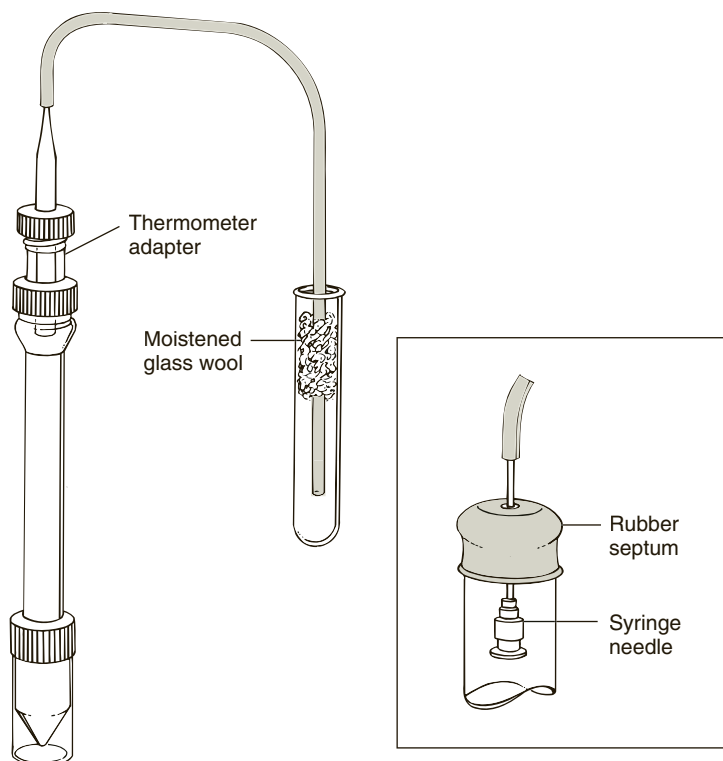


Figure 7.12

Microscale external gas trap. (The inset shows an expanded view of an alternative fitting, using a syringe needle and a rubber septum.)

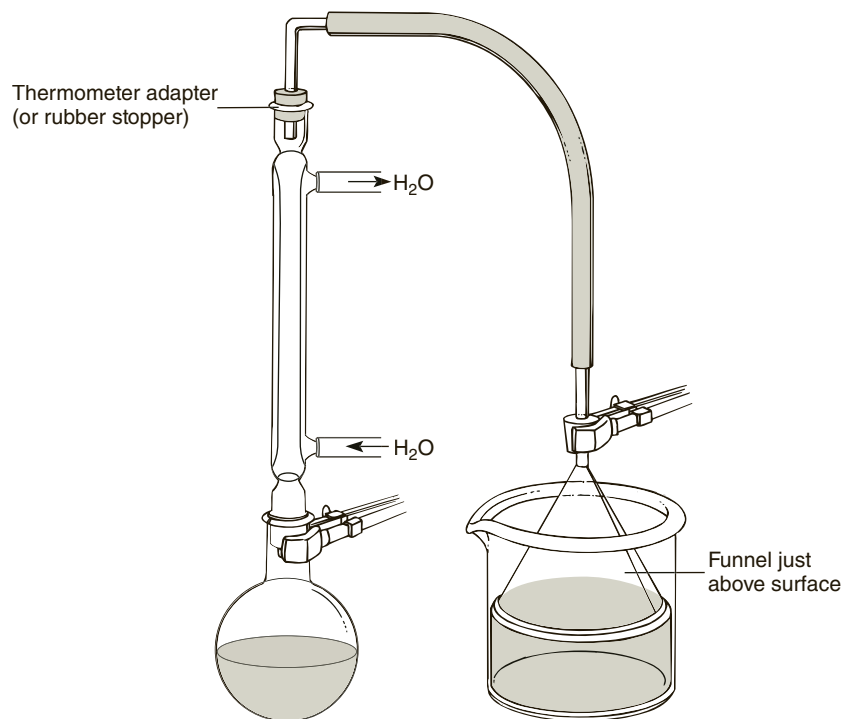


Figure 7.13
An inverted-funnel gas trap.

With large-scale reactions, a trap using an inverted funnel placed in a beaker of water is used. A piece of glass tubing, inserted through a thermometer adapter attached to the reaction apparatus, is connected to flexible tubing. The tubing is attached to a conical funnel. The funnel is clamped in place inverted over a beaker of water. The funnel is clamped so that its lip *almost touches* the water surface, but is not placed below the surface of the water. With this arrangement, water cannot be sucked back into the reaction if the pressure in the reaction vessel changes suddenly. This type of trap can also be used in microscale applications. An example of the inverted-funnel type of gas trap is shown in Figure 7.13.

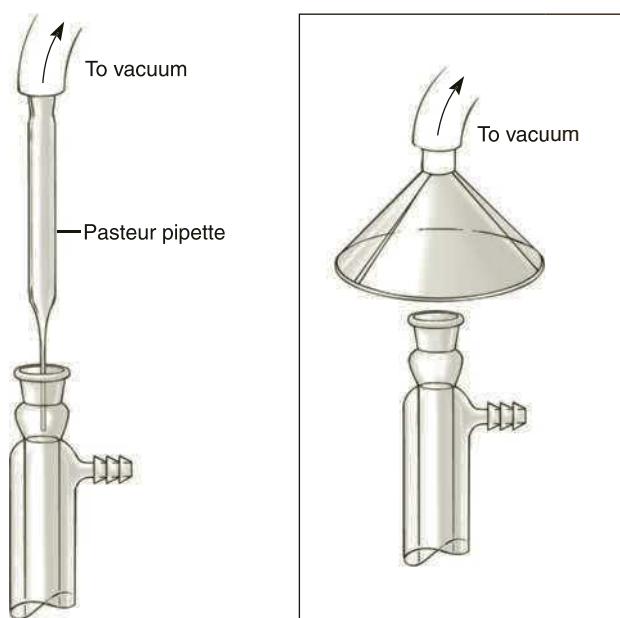


Figure 7.14
Removal of noxious gases under vacuum. (The inset shows an alternative assembly, using a funnel in place of the Pasteur pipette.)

C. Removal of Noxious Gases Using an Aspirator

An aspirator can be used to remove noxious gases from the reaction. The simplest approach is to clamp a disposable Pasteur pipette so that its tip is placed well into the condenser atop the reaction vial. An inverted funnel clamped over the apparatus can also be used. The pipette or funnel is attached to an aspirator with flexible tubing. A trap should be placed between the pipette or funnel and the aspirator. As gases are liberated from the reaction, they rise into the condenser. The vacuum draws the gases away from the apparatus. Both types of systems are shown in Figure 7.14. In the special case in which

the noxious gases are soluble in water, connecting a water aspirator to the pipette or funnel removes the gases from the reaction and traps them in the flowing water without the need for a separate gas trap.

7.9 Collecting Gaseous Products

In Section 7.8, means of removing unwanted gaseous products from the reaction system were examined. Some experiments produce gaseous products that you must collect and analyze. Methods to collect gaseous products are all based on the same principle. The gas is carried through tubing from the reaction to the opening of a flask or a test tube, which has been filled with water and is inverted in a container of water. The gas is allowed to bubble into the inverted collection tube (or flask). As the collection tube fills with gas, the water is displaced into the water container. If the collection tube is graduated, as in a graduated cylinder or a centrifuge tube, you can monitor the quantity of gas produced in the reaction.

If the inverted gas-collection tube is constructed from a piece of glass tubing, a rubber septum can be used to close the upper end of the container. This type of collection tube is shown in Figure 7.15. A sample of the gas can be removed using a syringe equipped with a needle. The gas that is removed can be analyzed by gas chromatography (see Technique 22).

Some of the glassware kits for microscale experiments contain a special, all-glass, capillary gas-delivery tube. The tube is attached to the top of the reaction apparatus by means of a ground-glass joint, and the open end of the capillary tubing is placed into an inverted, water-filled flask or test tube, clamped over a water bath. An

example of a microscale kit gas-delivery tube is shown in Figure 7.16A. This type of tube is an efficient means of collecting gases. A disadvantage, however, is that it is expensive and relatively easy to break.


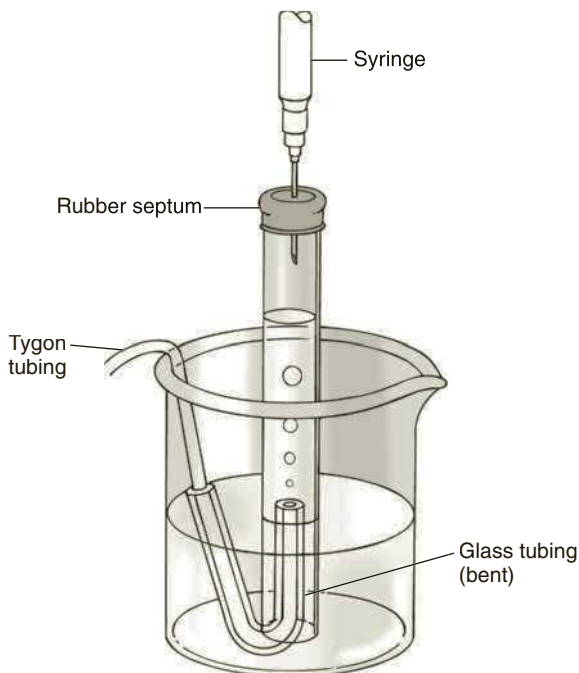


Figure 7.15
Gas-collection tube, with rubber septum.

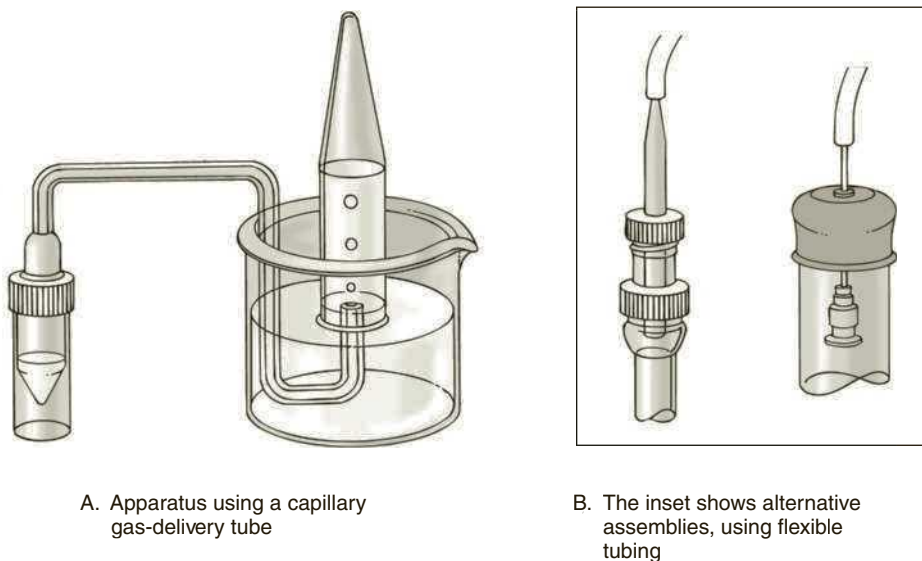

A. Apparatus using a capillary gas-delivery tube

B. The inset shows alternative assemblies, using flexible tubing

Figure 7.16
Gas-delivery tubes.

A simpler, less expensive approach is to use flexible tubing of a fine diameter to lead the gases from the reaction vessel to the collecting container. One method is to place a hypodermic syringe needle, point upward, through a rubber septum. The septum is attached to the top of the reaction apparatus, and a piece of fine flexible tubing is fitted over the end of the needle. The free end of the tubing is placed in the water bath, underneath the opening of the water-filled collection container. The gases bubble into the container, where they are collected. This alternative apparatus is shown in Figure 7.15 and also as an inset in Figure 7.16B.

Another alternative, which may also be used with larger-scale experiments, is to place a piece of glass tubing or the tip of a Pasteur pipette through a thermometer adapter. The thermometer adapter is attached to the top of the reaction apparatus, and flexible tubing is attached to the piece of glass tubing. The free end of the tubing is positioned in the opening of the water-filled collection vessel, as described previously. This variation is also shown as an inset in Figure 7.16. As an option, you may attach a second piece of glass tubing to the free end of the flexible hose. This piece of glass tubing sometimes makes it easier to fix the open end in the proper position in the opening of the collection flask.

7.10 Evaporation of Solvents

In many experiments, it is necessary to remove excess solvent from a solution. An obvious approach is to allow the container to stand unstoppered in the hood for several hours until the solvent has evaporated. This method is generally not practical, however, and a quicker, more efficient means of evaporating solvents must be used. Figures 7.17 and 7.18 show several methods of removing solvents by evaporation. Figure 7.17 depicts microscale methods; Figure 7.18 is devoted to large-scale procedures.

NOTE: It is good laboratory practice to evaporate solvents in the hood.

Microscale Methods. A simple means of evaporating a solvent is to place a conical vial in a warm water bath or a warm sand bath. The heat from the water bath or sand bath will warm the solvent to a temperature where it can evaporate within a short time. The heat from the water bath or sand bath can be adjusted to provide the best rate of evaporation, but the liquid should not be allowed to boil vigorously. The evaporation rate can be increased by allowing a stream of dry air or nitrogen to be directed into the vial (Figure 7.17A). The moving gas stream will sweep the vapors from the vial and accelerate the evaporation. As an alternative, a vacuum can be applied above the vial to draw away solvent vapors (Figure 7.17B and 7.17C).

A convenient water bath suitable for microscale methods can be constructed by placing the aluminum collars, which are generally used with aluminum heating blocks into a 150-mL beaker (Figure 7.17A). In some cases, it may be necessary to round off the sharp edges of the collars with a file in order to allow them to fit properly into the beaker. Held by the aluminum collars, the conical vial will stand securely in the beaker. This assembly can be filled with water and placed on a hot plate for use in the evaporation of small amounts of solvent.

Aluminum heating blocks placed on a hot plate can also be used for the evaporation of solvents (Figure 7.17B). You must be careful, however, not to allow the aluminum block to become too hot, or the sample may decompose thermally.

During a crystallization procedure, you often must remove excess solvent from the solution. If a Craig tube is being used for the crystallization, the excess solvent can be removed directly from the Craig tube (see Technique 11, Section 11.4). The Craig tube is placed in a warm water bath or warm sand bath. Alternatively,

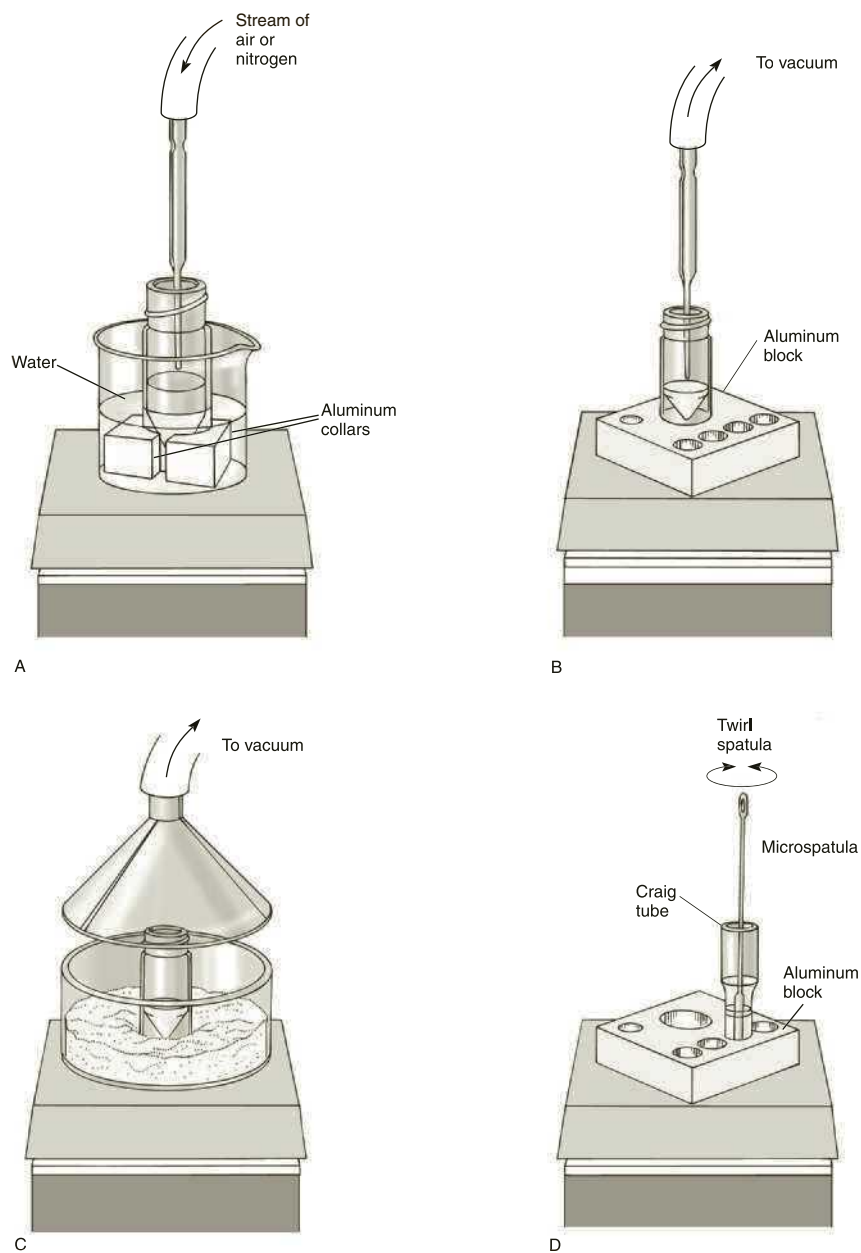


Figure 7.17
Evaporation of solvents (microscale methods).

the Craig tube can be placed into one of the small holes of an aluminum block. A microspatula is placed into the Craig tube, and it is twirled rapidly as the solvent evaporates (Figure 7.17D). The twirling spatula acts in the same manner as a boiling stone; it prevents bumping and accelerates the evaporation.

Commercially available evaporation stations may be useful when a large number of evaporations must be performed at the same time. This type of equipment consists of several holders for vials or flasks. At each position, a piece of tubing equipped with a metal tip is used to direct a stream of air into the vessel. A water bath is used to heat all the containers simultaneously.

Larger-Scale Methods. On a large scale, these evaporation methods can also be applied to standard-sized glassware. Solvents can be evaporated from solutions

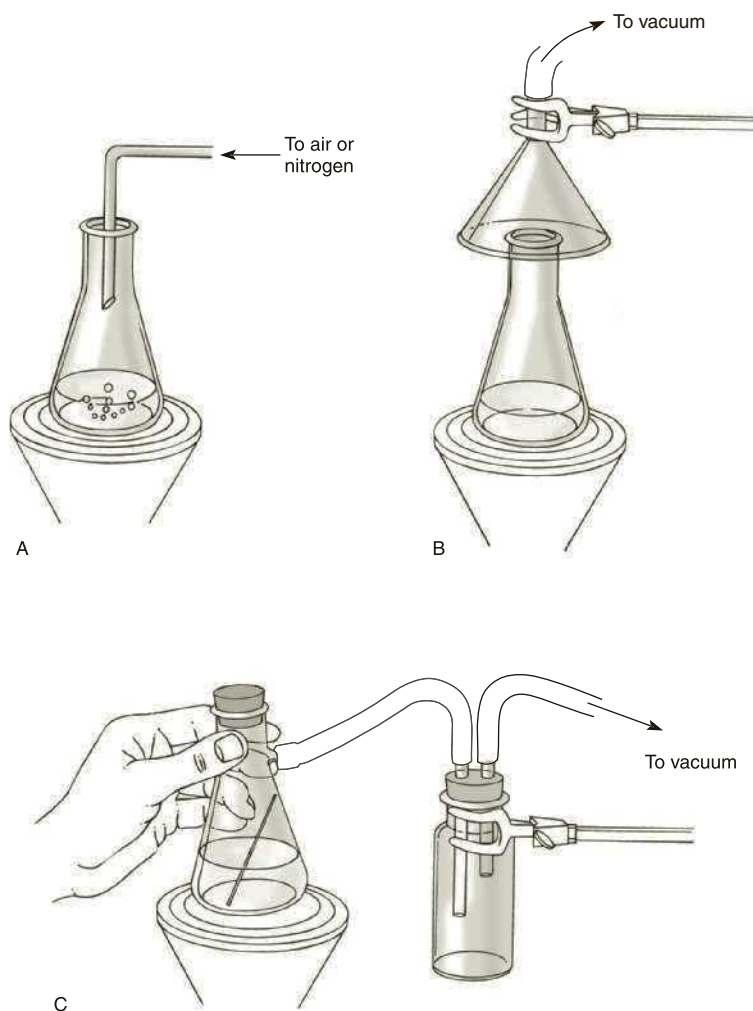


Figure 7.18
Evaporation of solvents (large-scale methods).

in Erlenmeyer flasks by adapting the techniques described previously. An Erlenmeyer flask can be placed on a source of heat, and the solvent can be removed by evaporation under a gas stream or a vacuum. Sources of heat that can be used with Erlenmeyer flasks include sand and steam baths and hot plates. A solution can also be placed in a side arm test tube or a filter flask, which is attached to a source of vacuum. A wooden stick or a piece of a melting point capillary is often placed in the solution, and the flask or test tube is swirled over the source of heat to reduce the possibility of bumping. The methods are illustrated in Figure 7.18.

7.11 Rotary Evaporator

In some organic chemistry laboratories, solvents are evaporated under reduced pressure using a **rotary evaporator**. This is a motor-driven device that is designed for rapid evaporation of solvents, with heating, while minimizing the possibility of bumping. A vacuum is applied to the flask, and the motor spins the flask. The rotation of the flask spreads a thin film of the liquid over the surface of the glass, which accelerates evaporation. The rotation also agitates the solution sufficiently to reduce the problem of bumping. A water bath can be placed under the flask to warm the solution and increase the vapor pressure of the solvent. One can select the speed at which the flask is rotated and the temperature of the water bath to attain the desired evaporation rate. As the solvent evaporates from the rotating flask,

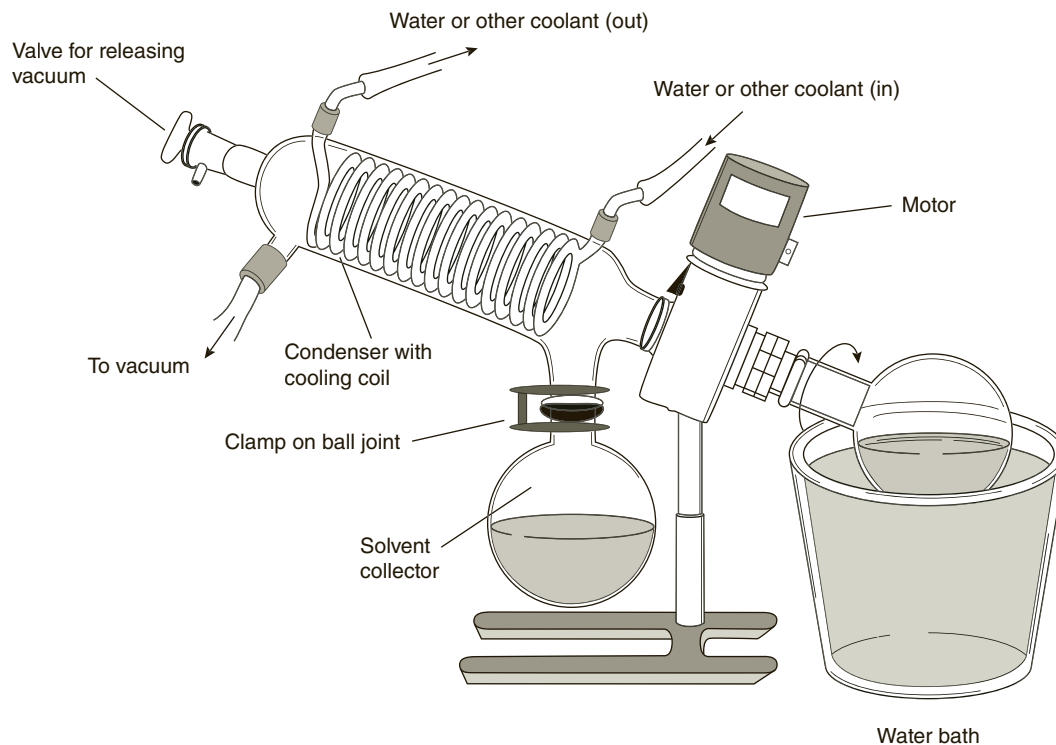


Figure 7.19
A rotary evaporator.

the vapors are cooled by the condenser, and the resulting liquid collects in the flask. The product remains behind in the rotating flask. A complete rotary evaporator assembly is shown in Figure 7.19. If the coolant is sufficiently cold, virtually all of the solvent can be recovered and recycled. This is a good example of *Green Chemistry* (see Green Chemistry essay, that precedes Experiment 28).

7.12 Microwave-Assisted Organic Chemistry

We are all familiar with the use of a microwave oven in the kitchen and its particular advantages. Cooking food in a microwave oven is much faster than in a conventional oven. Microwave cooking is much simpler, does not require as much crockery, and energy is not wasted in heating the container.

All of these advantages can also be applied to the chemistry laboratory. It is possible to conduct chemical reactions in much less time than with ordinary laboratory methods. Since the mid-1980s, chemists have been working on developing methods to apply microwave heating to chemical synthesis. Microwave-assisted organic chemical methods, or **microwave chemistry**, have gained wide acceptance, especially in industrial and research laboratories. Microwave heating is able to heat the chemical reagents without wasting energy in heating their container. In “green chemistry” applications, it allows the chemist to perform chemical reactions using less energy, in less time, often using water as a solvent, and often without using any solvent at all.

There does not seem to be general agreement as to the mechanism of microwave heating. The arguments are too complex to be included here. A basic understanding is possible, however. Microwave radiation is a form of *electromagnetic radiation*; this means that microwave radiation consists of oscillating electric and magnetic fields. When an oscillating electric field passes through a medium that contains polar or ionic substances, these molecules will attempt to orient themselves or oscillate in response to the electric field. Because these molecules are bound to surrounding



Figure 7.20
A microwave reaction system. (Reprinted courtesy of CEM Corporation.)

molecules in the medium, however, their motions are restricted, and they cannot respond completely to the oscillations of the electric field. This causes a non-equilibrium condition that results in an elevated instantaneous temperature in the immediate microscopic region surrounding the molecules that are being affected. As this localized temperature increases, molecules are activated above the required energy-of-activation threshold. Rates of reactions are dependent upon temperature; as the localized temperature increases, the molecules in that microscopic region will react faster.

Chemists first tried using domestic kitchen microwave ovens to speed up chemical reactions. They found that they were able to accelerate reactions, increase yields, and initiate otherwise impossible reactions. The results were often unsatisfactory, however, owing to uneven heating, lack of reproducibility, and the possibility

of explosions. The power output of a typical kitchen microwave oven cannot be adjusted. The oven cycles between periods of full power and periods of zero power. This means that the amount of microwave energy being transmitted into an experiment cannot be controlled precisely.


In recent years, companies have developed state-of-the-art microwave reaction systems to overcome these deficiencies. A modern reaction system, such as the one shown in Figure 7.20, has a specially designed vessel that focuses the microwave energy for efficient heating. Such systems are often equipped with automatic stirring and computer controls. Often a pressure control system may be included; this allows one to conduct a reaction at elevated temperature and pressure in the presence of volatile solvents or reagents. An automated sample changer is a useful accessory; this allows the chemist to conduct a series of repeated experiments without having to spend time watching the system.

Papers describing the advantages of microwave chemistry are appearing with increasing frequency in the chemical literature. Examples of experiments that can be conducted using microwave reaction systems include esterifications, condensation reactions, hydrogenations, cycloadditions, and even peptide syntheses. Besides offering a versatile method of chemical synthesis, microwave reaction systems also include the advantages that many of the reactions can be conducted in water, rather than in harmful organic solvents, or even in the complete absence of solvent. This capability makes microwave chemistry an important tool in “green chemistry.”





PROBLEMS

1. What is the best type of stirring device to use for stirring a reaction that takes place in the following type of glassware?
 - a. A conical vial
 - b. A 10-mL round-bottom flask
 - c. A 250-mL round-bottom flask

2. Should you use a drying tube for the following reaction? Explain.

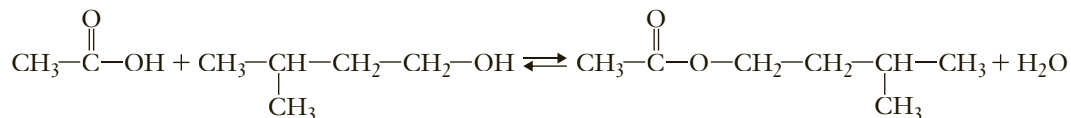

3. For which of the following reactions should you use a trap to collect noxious gases?

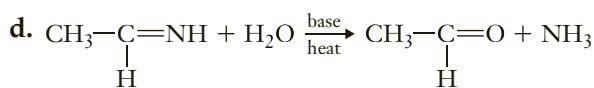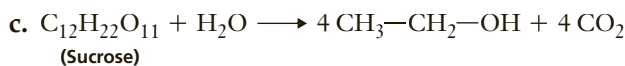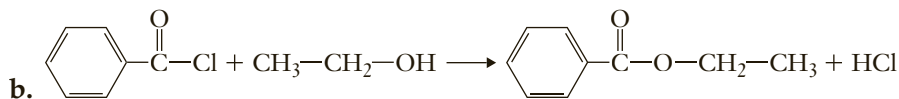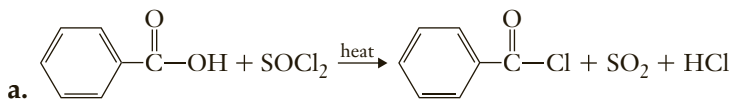

4. Criticize the following techniques:

- A reflux is conducted with a stopper in the top of the condenser.
- Water is passed through the reflux condenser at the rate of 1 gallon per minute.
- No water hoses are attached to the condenser during a reflux.
- A boiling stone is not added to the round-bottom flask until the mixture is boiling vigorously.
- To save money, you decide to save your boiling stones for another experiment.
- The reflux ring is located near the top of the condenser in a reflux setup.
- A rubber O-ring is omitted when the water condenser is attached to a conical vial.
- A gas trap is assembled with the funnel in Figure 7.13, completely submerged in the water in the beaker.
- Powdered drying agent is used rather than granular material.
- A reaction involving hydrogen chloride is conducted on the laboratory bench and not in a hood.
- An air-sensitive reaction apparatus is set up as shown in Figure 7.6.
- Air is used to evaporate solvent from an air-sensitive compound.

Filtration

Filtration is a technique used for two main purposes. The first is to remove solid impurities from a liquid. The second is to collect a desired solid from the solution from which it was precipitated or crystallized. Several different kinds of filtration are commonly used: two general methods include gravity filtration and vacuum (or suction) filtration. Two techniques specific to the microscale laboratory are filtration with a filter-tip pipette and filtration with a Craig tube. The various filtration techniques and their applications are summarized in Table 8.1. These techniques are discussed in more detail in the following sections.

Table 8.1 Filtration Methods

<i>Method</i>	<i>Application</i>	<i>Section</i>
Gravity filtration		
Filter cones	The volume of liquid to be filtered is about 10 mL or greater, and the solid collected in the filter is saved.	8.1A
Fluted filters	The volume of liquid to be filtered is greater than about 10 mL, and solid impurities are removed from a solution; often used in crystallization procedures.	8.1B
Filtering pipettes	Used with volumes less than about 10 mL to remove solid impurities from a liquid.	8.1C
Decantation	Although not a filtration technique, decantation can be used to separate a liquid from large, insoluble particles.	8.1D
Vacuum filtration		
Hirsch funnels	Used in the same way as Büchner funnels, except the volume of liquid is usually smaller (1–10 mL).	8.3
Büchner funnels	Primarily used to collect a desired solid from a liquid when the volume is greater than about 10 mL; used frequently to collect the crystals obtained from crystallization.	8.3
Filtering media	Used to remove finely divided impurities.	8.4
Filter-tip pipettes	May be used to remove a small amount of solid impurities from a small volume (1–2 mL) of liquid; also useful for pipetting volatile liquids, especially in extraction procedures.	8.6
Craig tubes	Used to collect a small amount of crystals resulting from crystallizations in which the volume of the solution is less than 2 mL.	8.7
Centrifugation	Although not strictly a filtration technique, centrifugation may be used to remove suspended impurities from a liquid (1–25 mL).	8.8

8.1 Gravity Filtration The most familiar filtration technique is probably filtration of a solution through a paper filter held in a funnel, allowing gravity to draw the liquid through the paper. Because even a small piece of filter paper will absorb a significant volume of liquid in most microscale procedures requiring filtration, this technique is useful only when the volume of mixture to be filtered is greater than 10 mL. For many macroscale and microscale procedures, a more suitable technique, which also makes use of gravity, is to use a Pasteur (or disposable) pipette with a cotton or glass wool plug (called a filtering pipette).

A. Filter Cones

This filtration technique is most useful when the solid material being filtered from a mixture is to be collected and used later. The filter cone, because of its smooth sides, can easily be scraped free of collected solids. Because of the many folds, fluted filter paper, described in the next section, cannot be scraped easily. The filter cone is likely to be used in microscale experiments only when a relatively large volume (greater than 10 mL) is being filtered and when a Hirsch funnel (Section 8.3) is not appropriate.

The filter cone is prepared as indicated in Figure 8.1. It is then placed into a funnel of an appropriate size. With filtrations using a simple filter cone, solvent may form seals between the filter and the funnel and between the funnel and the lip of the receiving flask. When a seal forms, the filtration stops because the displaced air has no possibility of escaping. To avoid the solvent seal, you can insert a small piece of paper, a paper clip, or some other bent wire between the funnel and the lip of the flask to let the displaced air escape. As an alternative, you can support the funnel by a clamp fixed *above* the flask rather than placed on the neck of the flask. A gravity filtration using a filter cone is shown in Figure 8.2.

B. Fluted Filters

This filtration method is also most useful when filtering a relatively large amount of liquid. Because a fluted filter is used when the desired material is expected to remain in solution, this filter is used to remove undesired solid materials, such as dirt particles, decolorizing charcoal, and undissolved impure crystals. A fluted filter is often used to filter a hot solution saturated with a solute during a crystallization procedure.

The technique for folding a fluted filter paper is shown in Figure 8.3. An advantage of a fluted filter is that it increases the speed of filtration in two ways. First, it increases the surface area of the filter paper through which the solvent

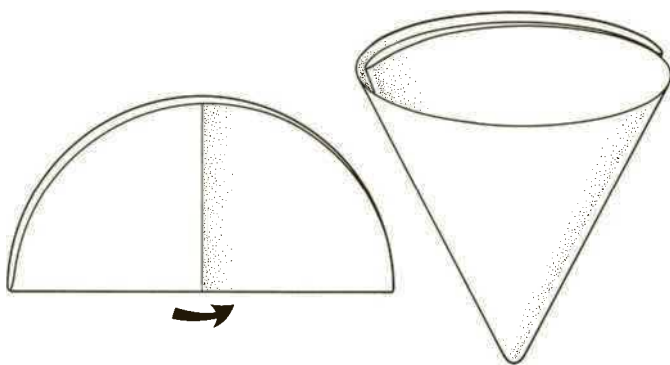


Figure 8.1
Folding a filter cone.

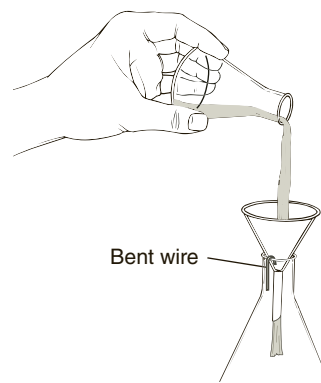
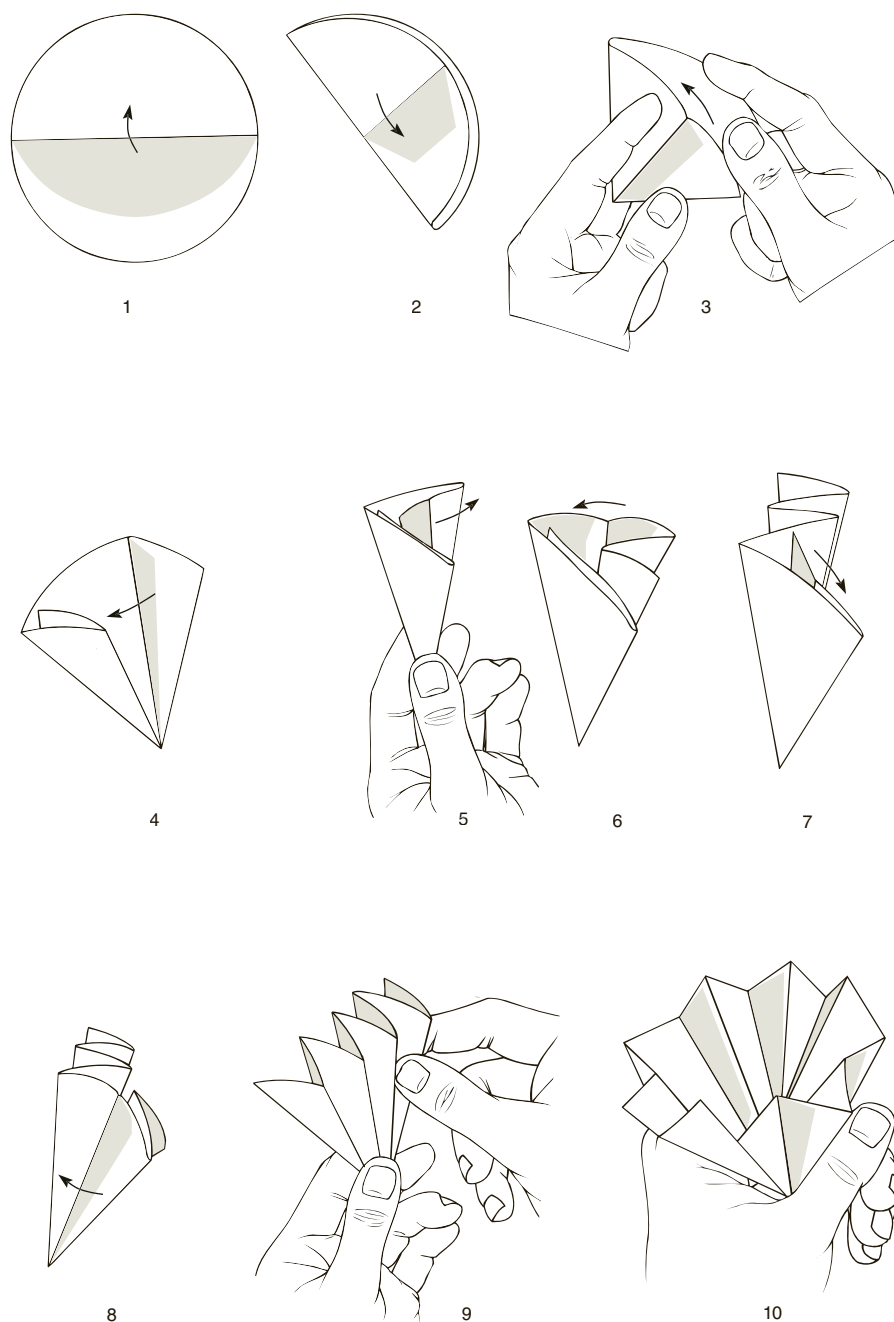


Figure 8.2
Gravity filtration with a filter cone.

**Figure 8.3**

Folding a fluted filter paper, or origami at work in the organic chemistry laboratory.

seeps; second, it allows air to enter the flask along its sides to permit rapid pressure equalization. If pressure builds up in the flask from hot vapors, filtering slows down. This problem is especially pronounced with filter cones. The fluted filter tends to reduce this problem considerably, but it may be a good idea to clamp the funnel above the receiving flask or to use a piece of paper, paper clip, or wire between the funnel and the lip of the flask as an added precaution against solvent seals.

Filtration with a fluted filter is relatively easy to perform when the mixture is at room temperature. However, when it is necessary to filter a hot solution saturated

with a dissolved solute, a number of steps must be taken to ensure that the filter does not become clogged by solid material accumulated in the stem of the funnel or in the filter paper. When the hot, saturated solution comes in contact with a relatively cold funnel (or a cold flask, for that matter), the solution is cooled and may become supersaturated. If crystallization then occurs in the filter, either the crystals will fail to pass through the filter paper or they will clog the stem of the funnel.

To keep the filter from clogging, use one of the following four methods. The first is to use a short-stemmed or a stemless funnel. With these funnels, it is less likely that the stem of the funnel will become clogged by solid material. The second method is to keep the liquid to be filtered at or near its boiling point at all times. The third way is to preheat the funnel by pouring hot solvent through it before the actual filtration. This keeps the cold glass from causing instantaneous crystallization. And fourth, it is helpful to keep the **filtrate** (filtered solution) in the receiver hot enough to continue boiling *slightly* (by setting it on a hot plate, for example). The refluxing solvent heats the receiving flask and the funnel stem and washes them clean of solids. This boiling of the filtrate also keeps the liquid in the funnel warm.

C. Filtering Pipettes

A filtering pipette is a microscale technique most often used to remove solid impurities from a liquid with a volume less than 10 mL. It is important that the mixture being filtered be at or near room temperature because it is difficult to prevent premature crystallization in a hot solution saturated with a solute.

To prepare this filtration device, a small piece of cotton is inserted into the top of a Pasteur (disposable) pipette and pushed down to the beginning of the lower constriction in the pipette, as shown in Figure 8.4. It is important to use enough cotton to collect all the solid being filtered; however, the amount of cotton used should not be so large that the flow rate through the pipette is significantly restricted. For the

same reason, the cotton should not be packed too tightly. The cotton plug can be pushed down gently with a long thin object such as a glass stirring rod or a wooden applicator stick. It is advisable to wash the cotton plug by passing about 1 mL of solvent (usually the same solvent that is to be filtered) through the filter.

In some cases, such as when filtering a strongly acidic mixture or when performing a very rapid filtration to remove dirt or impurities of large particle size from a solution, it may be better to use glass wool in place of the cotton. The disadvantage in using glass wool is that the fibers do not pack together as tightly, and small particles will pass through the filter more easily.

To conduct a filtration (with either a cotton or glass wool plug), the filtering pipette is clamped so that the filtrate will drain into an appropriate container. The mixture to be filtered is usually transferred to the filtering pipette with another Pasteur pipette. If a small volume of liquid is being filtered (less than 1 mL or 2 mL), it is advisable to rinse the filter and plug with a small amount of solvent after the last of the filtrate has passed through the filter. The rinse solvent is then combined with the

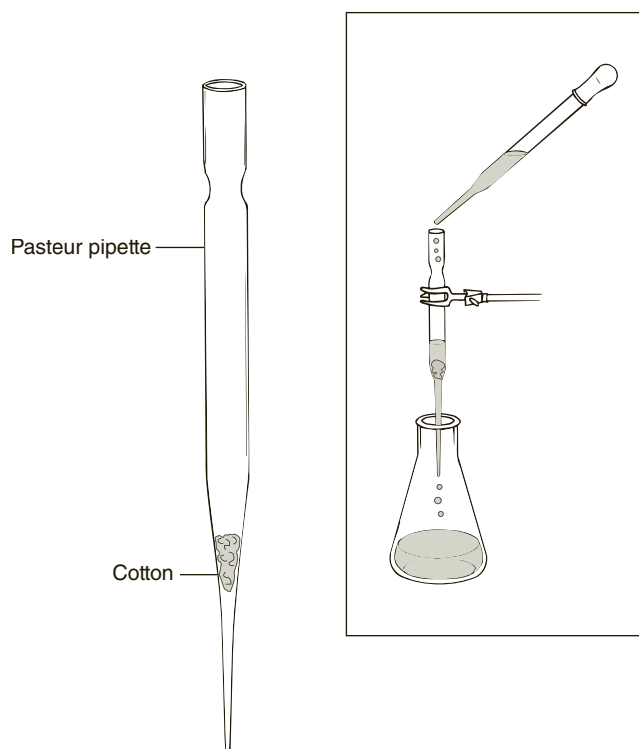


Figure 8.4
A filtering pipette.

original filtrate. The rate of filtration can be increased by gently applying pressure to the top of the pipette using a pipette bulb, if desired.

Depending on the amount of solid being filtered and the size of the particles (small particles are more difficult to remove by filtration), it may be necessary to put the filtrate through a second filtering pipette. This should be done with a new filtering pipette rather than with the one already used.

D. Decantation

It is not always necessary to use filter paper to separate insoluble particles. If you have large, heavy, insoluble particles, with careful pouring you can decant the solution, leaving behind the solid particles that will settle to the bottom of the flask. The term *decant* means “to carefully pour out the liquid, leaving the insoluble particles behind.” For example, boiling stones or sand granules in the bottom of an Erlenmeyer flask filled with a liquid can easily be separated in this way. This procedure is often preferred over filtration and usually results in a smaller loss of material. If there are a large number of particles and they retain a significant amount of the liquid, they can be rinsed with solvent and a second decantation performed. The term *decant* was coined in the wine industry, where it is often necessary to let the wine settle and then carefully pour it out of the original bottle into a clean one, leaving the “must” (insoluble particles) behind.

8.2 Filter Paper

Many kinds and grades of filter paper are available. The paper must be correct for a given application. In choosing filter paper, you should be aware of its various properties. **Porosity** is a measure of the size of the particles that can pass through the paper. Highly porous paper does not remove small particles from solution; paper with low porosity removes very small particles. **Retentivity** is a property that is the opposite of porosity. Paper with low retentivity does not remove small particles from the filtrate. The **speed** of filter paper is a measure of the time it takes a liquid to drain through the filter. Fast paper allows the liquid to drain quickly; with slow paper, it takes much longer to complete the filtration. Because all these properties are related, fast filter paper usually has a low retentivity and high porosity, and slow filter paper usually has high retentivity and low porosity.

Table 8.2 compares some commonly available qualitative filter paper types and ranks them according to porosity, retentivity, and speed. Eaton–Dikeman (E&D),

Table 8.2 Some common qualitative filter paper types and approximate relative speeds and retentivities

	Fine	High	Slow	Type (by Number)			
				E&D	S&S	Whatman	
Porosity	↓	↓	↓	Speed			
				Very slow	610	576	5
				Slow	613	602	3
				Medium	615	597	2
				Fast	617	595	1
				Very fast	—	604	4
	Coarse	Low	Fast				

Schleicher and Schuell (S&S), and Whatman are the most common brands of filter paper. The numbers in the table refer to the grades of paper used by each company.

8.3 Vacuum Filtration

Vacuum, or suction, filtration is more rapid than gravity filtration and is most often used to collect solid products resulting from precipitation or crystallization. This technique is used primarily when the volume of liquid being filtered is more than 1–2 mL. With smaller volumes, use of the Craig tube (Section 8.7) is the preferred technique. In a vacuum filtration, a receiver flask with a side arm, a **filter flask**, is used. For microscale work, the most useful size is a 50-mL filter flask. For macroscale laboratory work, the most useful sizes of filter flasks range from 50 mL to 500 mL, depending on the volume of liquid being filtered. The side arm is connected by *heavywalled* rubber tubing (see Technique 16, Figure 16.3) to a source of vacuum. Thin-walled tubing will collapse under vacuum, due to atmospheric pressure on its outside walls, and will seal the vacuum source from the flask. Because this apparatus is unstable and can tip over easily, it must be clamped, as shown in Figure 8.5.

CAUTION



It is essential that the filter flask be clamped.

Two types of funnels are useful for vacuum filtration, the Hirsch funnel and the Büchner funnel. The **Hirsch funnel** is used for filtering smaller amounts of solid from solution. Hirsch funnels are usually made from polypropylene or porcelain. The polypropylene Hirsch funnel (see Figure 8.5A) is sealed to a 50-mL filter flask by a small section of Gooch tubing. This Hirsch funnel has a built-in adapter that forms a tight seal with some 25-mL filter flasks without the Gooch tubing. A fritted polyethylene

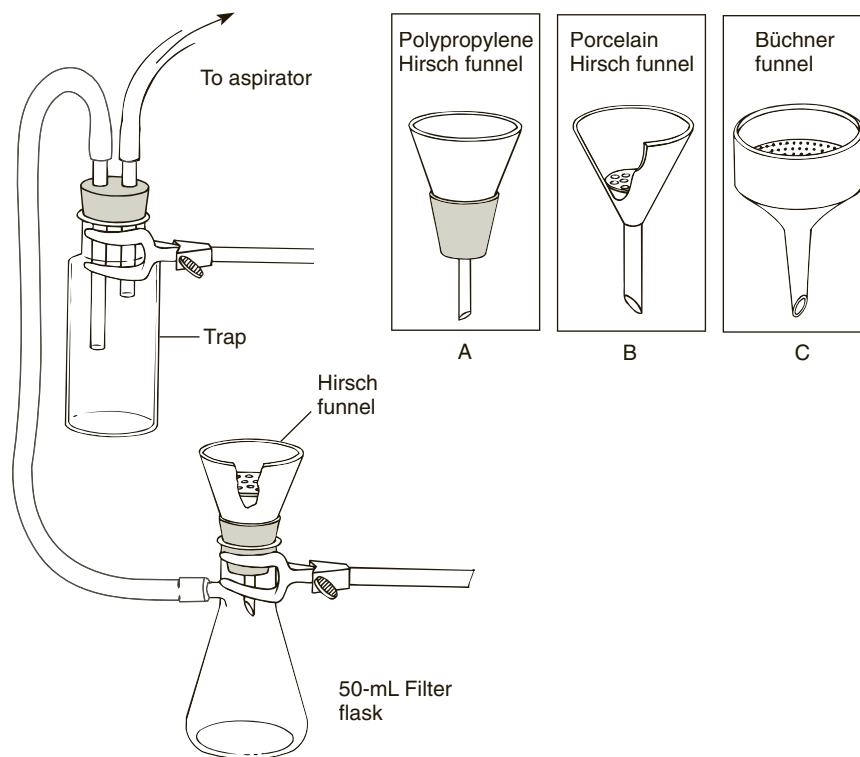


Figure 8.5
Vacuum filtration.

disk fits into the bottom of the funnel. To prevent the holes in this disk from becoming clogged with solid material, the funnel should always be used with a circular filter paper that has the same diameter (1.27 cm) as the polyethylene disk. Before beginning the filtration, it is advisable to moisten the paper with a small amount of solvent. The moistened filter paper adheres more strongly to the fritted disk and prevents unfiltered mixture from passing around the edges of the filter paper. A porcelain Hirsch funnel (see Figure 8.5B) is sealed to the filter flask by a rubber stopper or a filter (Neoprene) adapter. The flat bottom of this Hirsch funnel, which should be 1–2 cm in diameter, is covered with an unfolded piece of circular filter paper. To prevent the escape of solid materials from the funnel, you must be certain that the filter paper fits the funnel exactly. It must cover all the holes in the bottom of the funnel but not extend up the sides of the funnel. With a porcelain Hirsch funnel, it is also important to moisten the paper with a small amount of solvent before beginning the filtration.

The **Büchner funnel**, which is shown in Figure 8.5C, operates on the same principle as the Hirsch funnel, but it is usually larger, and its sides are vertical rather than sloped. It is sealed to the filter flask with a rubber stopper or a Neoprene adapter. In the Büchner funnel, the filter paper must also cover all the holes in the bottom but must not extend up the sides.

Because the filter flask is attached to a source of vacuum, a solution poured into a Hirsch funnel or Büchner funnel is literally sucked rapidly through the filter paper. For this reason, vacuum filtration is generally not used to separate fine particles such as decolorizing charcoal, because the small particles would likely be pulled through the filter paper. However, this problem can be alleviated, when desired, by the use of specially prepared filter beds (see Section 8.4).

8.4 Filtering Media

It is occasionally necessary to use specially prepared filter beds to separate fine particles when using vacuum filtration. Often, very fine particles either pass right through a paper filter or clog it so completely that the filtering stops. This is avoided by using a substance called Filter Aid, or Celite. This material is also called **diatomaceous earth** because of its source. It is a finely divided inert material derived from the microscopic shells of dead diatoms (a type of phytoplankton that grows in the sea).

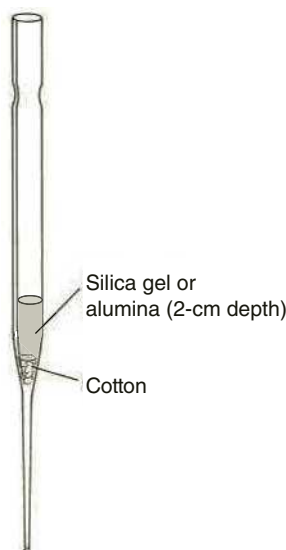


Figure 8.6
A Pasteur pipette with filtering media.

CAUTION



Diatomaceous earth is a lung irritant. When using Filter Aid, take care not to breathe the dust.

Filter Aid will not clog the fiber pores of filter paper. It is **slurried**, mixed with a solvent to form a rather thin paste, and filtered through a Hirsch or Büchner funnel (with filter paper in place) until a layer of diatoms about 2–3 mm thick is formed on top of the filter paper. The solvent in which the diatoms were slurried is poured from the filter flask, and, if necessary, the filter flask is cleaned before the actual filtration is begun. Finely divided particles can now be suction-filtered through this layer and will be caught in the Filter Aid. This technique is used for removing impurities, not for collecting a product. The filtrate (filtered solution) is the desired material in this procedure. If the material caught in the filter were the desired material, you would have to try to separate the product from all those diatoms! Filtration with Filter Aid is not appropriate when the desired substance is likely to precipitate or crystallize from solution.

In microscale work, it may sometimes be more convenient to use a column prepared with a Pasteur pipette to separate fine particles from a solution. The Pasteur pipette is packed with alumina or silica gel, as shown in Figure 8.6.

8.5 The Aspirator

The most common source of vacuum (approximately 10–20 mm Hg) in the laboratory is the water aspirator, or “water pump,” illustrated in Figure 8.7. This device passes water rapidly past a small hole to which a side arm is attached. The water pulls air in through the side arm. This phenomenon, called the Bernoulli effect, causes a reduced pressure along the side of the rapidly moving water stream and creates a partial vacuum in the side arm.

NOTE: The aspirator works most effectively when the water is turned on to the fullest extent.

A water aspirator can never lower the pressure beyond the vapor pressure of the water used to create the vacuum. Hence, there is a lower limit to the pressure (on cold days) of 9–10 mm Hg. A water aspirator does not provide as high a vacuum in the summer as in the winter, because of this water-temperature effect.

A trap must be used with an aspirator. One type of trap is illustrated in Figure 8.5. Another method for securing this type of trap is shown in Figure 8.8. This simple holder can be constructed from readily available material and can be placed anywhere on the laboratory bench. Although not often needed, a trap can prevent water from contaminating your experiment. If the water pressure in the laboratory drops suddenly, the pressure in the filter flask may suddenly become lower than the pressure in the water aspirator. This would cause water to be drawn from the aspirator stream into the filter flask and contaminate the filtrate or even the material in the filter. The trap stops this reverse flow. A similar flow will occur if the water flow at the aspirator is stopped before the tubing connected to the aspirator side arm is disconnected.

NOTE: Always disconnect the tubing before stopping the aspirator.

If a “backup” begins, disconnect the tubing as rapidly as possible before the trap fills with water. Some chemists like to fit a stopcock into the stopper on top of the trap. A three-hole stopper is required for this purpose. With a stopcock in the trap, the system can be vented before the aspirator is shut off. Then water cannot back up into the trap.

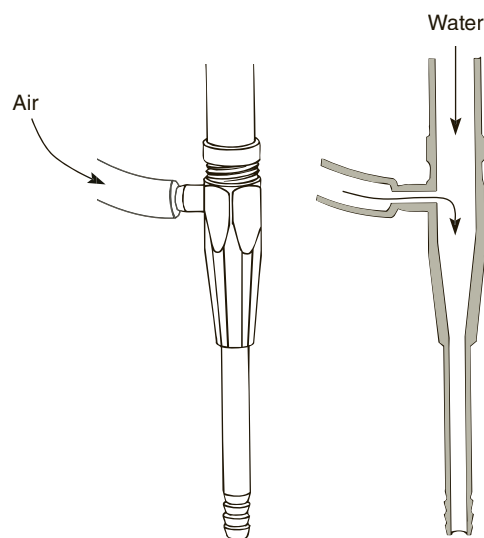


Figure 8.7
An aspirator.

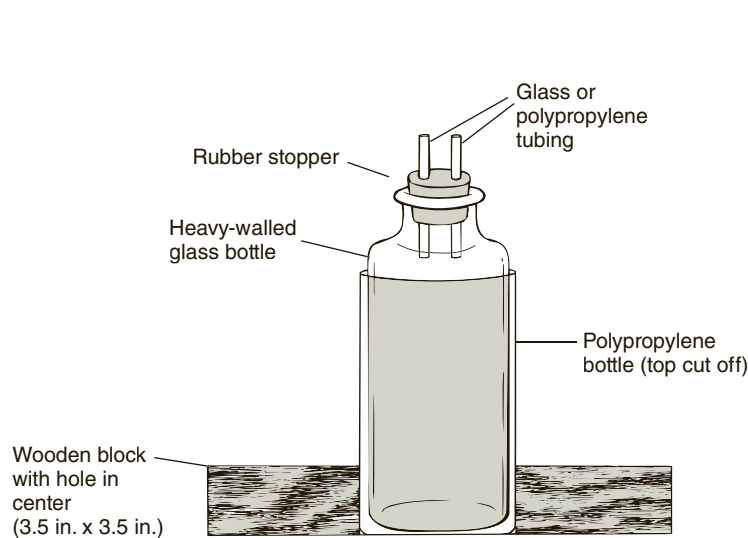


Figure 8.8
A simple aspirator trap and holder.

Aspirators do not work well if too many people use the water line at the same time because the water pressure is lowered. Also, the sinks at the ends of the lab benches or the lines that carry away the water flow may have a limited capacity for draining the resultant water flow from too many aspirators. Care must be taken to avoid floods.

8.6 Filter-Tip Pipette

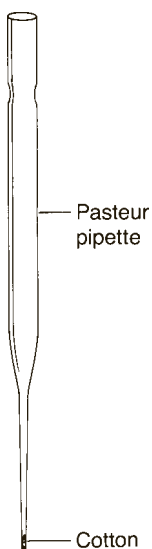


Figure 8.9
A filter-tip pipette.

The filter-tip pipette, illustrated in Figure 8.9, has two common uses. The first is to remove a small amount of solid, such as dirt or filter paper fibers, from a small volume of liquid (1–2 mL). It can also be helpful when using a Pasteur pipette to transfer a highly volatile liquid, especially during an extraction procedure (see Technique 12, Section 12.5).

Preparing a filter-tip pipette is similar to preparing a filtering pipette, except that a much smaller amount of cotton is used. A *tiny* piece of cotton is loosely shaped into a ball and placed into the large end of a Pasteur pipette. Using a wire with a diameter slightly smaller than the inside diameter of the narrow end of the pipette, push the ball of cotton to the bottom of the pipette. If it becomes difficult to push the cotton, you have probably started with too much cotton; if the cotton slides through the narrow end with little resistance, you probably have not used enough.

To use a filter-tip pipette as a filter, the mixture is drawn up into the Pasteur pipette using a pipette bulb and then expelled. With this procedure, a small amount of solid will be captured by the cotton. However, very fine particles, such as activated charcoal, cannot be removed efficiently with a filter-tip pipette, and this technique is not effective in removing more than a trace amount of solid from a liquid.

Transferring many organic liquids with a Pasteur pipette can be a somewhat difficult procedure for two reasons. First, the liquid may not adhere well to the glass. Second, as you handle the Pasteur pipette, the temperature of the liquid in the pipette increases slightly, and the increased vapor pressure may tend to “squirt” the liquid out the end of the pipette. This problem can be particularly troublesome when separating two liquids during an extraction procedure. The purpose of the cotton plug in this situation is to slow the rate of flow through the end of the pipette so you can control the movement of liquid in the Pasteur pipette more easily.

8.7 Craig Tubes

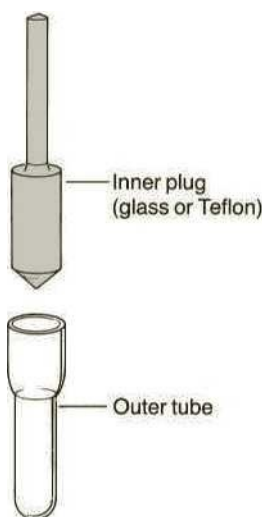


Figure 8.10
A Craig tube (2 mL).

The **Craig tube**, illustrated in Figure 8.10, is used primarily to separate crystals from a solution after a microscale crystallization procedure has been performed (Technique 11, Section 11.4). Although it may not be a filtration procedure in the traditional sense, the outcome is similar. The outer part of the Craig tube is similar to a test tube, except that the diameter of the tube becomes wider part of the way up the tube, and the glass is ground at this point so that the inside surface is rough. The inner part (plug) of the Craig tube may be made of Teflon or glass. If this part is glass, the end of the plug is also ground. With either a glass or a Teflon inner plug, there is only a partial seal where the plug and the outer tube come together. Liquid may pass through, but solid will not. This is the place where the solution is separated from the crystals.

After crystallization has been completed in the outer Craig tube, replace the inner plug (if necessary) and connect a thin copper wire or strong thread to the narrow part of the inner plug, as indicated in Figure 8.11A. While holding the Craig tube in an upright position, place a plastic centrifuge tube over the Craig tube so that the bottom of the centrifuge tube rests on top of the inner plug, as shown in Figure 8.11B. The copper wire should extend just below the lip of the centrifuge tube and is now bent upward around the lip of the centrifuge tube. This apparatus is then turned over so that the centrifuge tube is upright. The Craig tube is spun in a centrifuge (be sure it is balanced by placing another tube filled with water on the

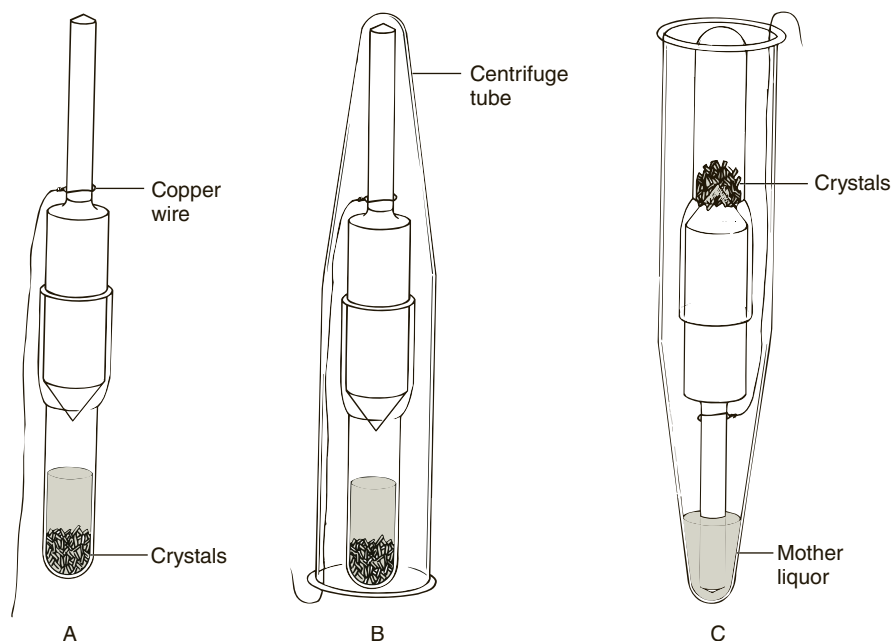


Figure 8.11
Separation with a Craig tube.

opposite side of the centrifuge) for several minutes until the **mother liquor** (solution from which the crystals grew) goes to the bottom of the centrifuge tube and the crystals collect on the end of the inner plug (see Figure 8.11C). Depending on the consistency of the crystals and the speed of the centrifuge, the crystals may spin down to the inner plug, or (if you are unlucky) they may remain at the other end of the Craig tube.¹ If the latter situation occurs, it may be helpful to centrifuge the Craig tube longer or, if this problem is anticipated, to stir the crystal- and solution-mixture with a spatula or stirring rod before centrifugation.

Using the copper wire, then pull the Craig tube out of the centrifuge tube. If the crystals collected on the end of the inner plug, it is now a simple procedure to remove the plug and scrape the crystals with a spatula onto a watch glass, a clay plate, or a piece of smooth paper. Otherwise, it will be necessary to scrape the crystals from the inside surface of the outer part of the Craig tube.

8.8 Centrifugation

Sometimes, centrifugation is more effective in removing solid impurities than conventional filtration techniques. Centrifugation is particularly effective in removing suspended particles, which are so small that the particles would pass through most filtering devices. Centrifugation may also be useful when the mixture must be kept hot to prevent premature crystallization while the solid impurities are removed.

Centrifugation is performed by placing the mixture in one or two centrifuge tubes (be sure to balance the centrifuge) and centrifuging for several minutes. The supernatant liquid is then decanted (poured off) or removed with a Pasteur pipette.

¹ *Note to the Instructor:* In some centrifuges, the bottom of the Craig tube may be close to the center of the centrifuge when the Craig tube assembly is placed into the centrifuge. In this situation, very little centrifugal force will be applied to the crystals, and it is likely that the crystals will not spin down. It may then be helpful to use an inner plug with a shorter stem. The stem on a Teflon inner plug can be easily cut off about 0.5 inch with a pair of wire cutters. This will help to spin down the crystals to the inner plug, and also the centrifuge can be run at a lower speed, which can help prevent breakage of the Craig tube.

PROBLEMS

1. In each of the following situations, what type of filtration device would you use? (Do not consider centrifugation when answering this question.)
 - a. Remove powdered decolorizing charcoal from 20 mL of solution
 - b. Collect crystals obtained from crystallizing a substance from about 1 mL of solution
 - c. Remove a very small amount of dirt from 1 mL of liquid
 - d. Isolate 0.2 g of crystals from about 5 mL of solution after performing a crystallization
 - e. Remove dissolved colored impurities from about 3 mL of solution
 - f. Remove solid impurities from 5 mL of liquid at room temperature

Physical Constants of Solids: The Melting Point

9.1 Physical Properties

The physical properties of a compound are those properties that are intrinsic to a given compound when it is pure. A compound may often be identified simply by determining a number of its physical properties. The most commonly recognized physical properties of a compound include its color, melting point, boiling point, density, refractive index, molecular weight, and optical rotation. Modern chemists would include the various types of spectra (infrared, nuclear magnetic resonance, mass, and ultraviolet-visible) among the physical properties of a compound. A compound's spectra do not vary from one pure sample to another. Here, we look at methods of determining the melting point. Boiling point and density of compounds are covered in Technique 13. Refractive index, optical rotation, and spectra are also considered separately.

Many reference books list the physical properties of substances. You should consult Technique 4 for a complete discussion on how to find data for specific compounds. The works most useful for finding lists of values for the nonspectroscopic physical properties include

The Merck Index

The CRC Handbook of Chemistry and Physics

Lange's Handbook of Chemistry

Aldrich Handbook of Fine Chemicals

Complete citations for these references can be found in Technique 29. Although the *CRC Handbook* has good tables, it adheres strictly to IUPAC nomenclature. For this reason, it may be easier to use one of the other references, particularly *The Merck Index* or the *Aldrich Handbook of Fine Chemicals*, in your first attempt to locate information (see Technique 4).

9.2 The Melting Point

The melting point of a compound is used by the organic chemist not only to identify the compound but also to establish its purity. A small amount of material is heated *slowly* in a special apparatus equipped with a thermometer or thermocouple, a heating bath or heating coil, and a magnifying eyepiece for observing the sample. Two temperatures are noted. The first is the point at which the first drop of liquid forms among the crystals; the second is the point at which the whole mass of crystals turns to a *clear* liquid. The melting point is recorded by giving this range of melting. You might say, for example, that the melting point of a substance is 51–54°C. That is, the substance melted over a 3-degree range.

The melting point indicates purity in two ways. First, the purer the material, the higher its melting point. Second, the purer the material, the narrower its melting-point range. Adding successive amounts of an impurity to a pure substance generally causes its melting point to decrease in proportion to the amount of

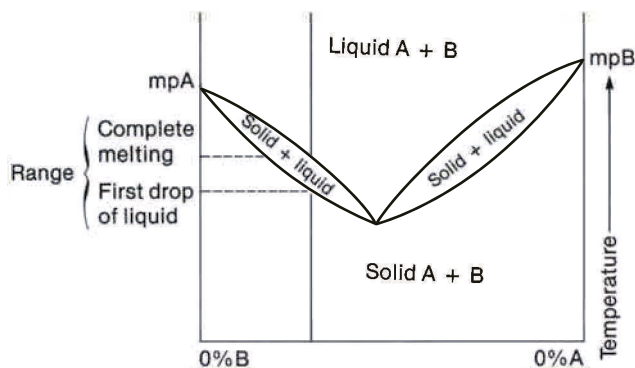


Figure 9.1
A melting-point-composition curve.

begin with pure A, the melting point decreases as impurity B is added. At some point, a minimum temperature, or **eutectic**, is reached, and the melting point begins to increase to that of substance B. The vertical distance between the lower and upper curves represents the melting range. Notice that for mixtures that contain relatively small amounts of impurity (<15%) and are not close to the eutectic, the melting range increases as the sample becomes less pure. The range indicated by the lines in Figure 9.1 represents the typical behavior.

We can generalize the behavior shown in Figure 9.1. Pure substances melt with a narrow range of melting. With impure substances, the melting range becomes wider, and the entire melting range is lowered. Be careful to note, however, that at the minimum point of the melting-point-composition curves, the mixture often forms a eutectic, which also melts sharply. Not all binary mixtures form eutectics, and some caution must be exercised in assuming that every binary mixture follows the previously described behavior. Some mixtures may form more than one eutectic; others might not form even one. In spite of these variations, both the melting point and its range are useful indications of purity, and they are easily determined by simple experimental methods.

9.3 Melting-Point Theory

Figure 9.2 is a phase diagram describing the usual behavior of a two-component mixture (A + B) on melting. The behavior on melting depends on the relative amounts of A and B in the mixture. If A is a pure substance (no B), then A melts sharply at its melting point t_A . This is represented by point A on the left side of the diagram. When B is a pure substance, it melts at t_B ; its melting point is represented by point B on the right side of the diagram. At either point A or point B, the pure solid passes cleanly, with a narrow range, from solid to liquid.

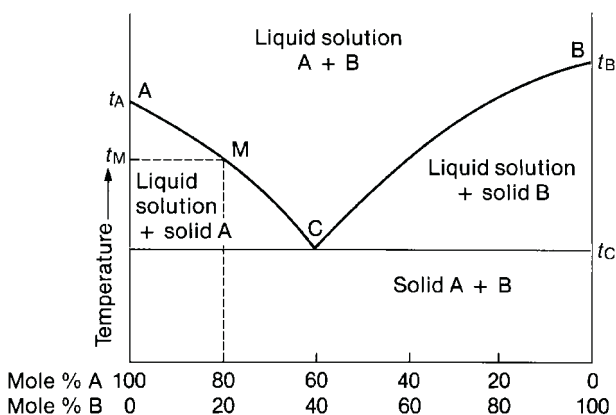


Figure 9.2
A phase diagram for melting in a two-component system.

In mixtures of A and B, the behavior is different. Using Figure 9.2, consider a mixture of 80% A and 20% B on a mole-per-mole basis (that is, mole percentage). The melting point of this mixture is given by t_M at point M on the diagram. That is, adding B to A has lowered the melting point of A from t_A to t_M . It has also expanded the melting range. The temperature t_M corresponds to the **upper limit** of the melting range.

Lowering the melting point of A by adding impurity B comes about in the following way. Substance A has the lower melting point in the phase diagram shown, and if heated, it begins to melt first. As A begins to melt, solid B begins to dissolve in the liquid A that is formed. When solid B dissolves in liquid A, the melting point is depressed. To understand this, consider the melting point from the opposite direction. When a liquid at a high temperature cools, it reaches a point at which it solidifies, or “freezes.” The temperature at which a liquid freezes is identical to its melting point. Recall that the freezing point of a liquid can be lowered by adding an impurity. Because the freezing point and the melting point are identical, lowering the freezing point corresponds to lowering the melting point. Therefore, as more impurity is added to a solid, its melting point becomes lower. There is, however, a limit to how far the melting point can be depressed. You cannot dissolve an infinite amount of the impurity substance in the liquid. At some point, the liquid will become saturated with the impurity substance. The solubility of B in A has an upper limit. In Figure 9.2, the solubility limit of B in liquid A is reached at point C, the **eutectic point**. The melting point of the mixture cannot be lowered below t_C , the melting temperature of the eutectic.

Now consider what happens when the melting point of a mixture of 80% A and 20% B is approached. As the temperature is increased, A begins to “melt.” This is not really a visible phenomenon in the beginning stages; it happens before liquid is visible. It is a softening of the compound to a point at which it can begin to mix with the impurity. As A begins to soften, it dissolves B. As it dissolves B, the melting point is lowered. The lowering continues until all B is dissolved or until the eutectic composition (saturation) is reached. When the maximum possible amount of B has been dissolved, actual melting begins, and one can observe the first appearance of liquid. The initial temperature of melting will be below t_A . The amount below t_A at which melting begins is determined by the amount of B dissolved in A but will never be below t_C . Once all B has been dissolved, the melting point of the mixture begins to rise as more A begins to melt. As more A melts, the semisolid solution is diluted by more A, and its melting point rises. While all this is happening, you can observe *both* solid and liquid in the melting-point capillary. Once all A has begun to melt, the composition of the mixture M becomes uniform and will reach 80% A and 20% B. At this point, the mixture finally melts sharply, giving a clear solution. The maximum melting-point range will be $t_C - t_M$, because t_A is depressed by the impurity B that is present. The lower end of the melting range will always be t_C ; however, melting will not always be observed at this temperature. An observable melting at t_C comes about only when a large amount of B is present. Otherwise, the amount of liquid formed at t_C will be too small to observe. Therefore, the melting behavior that is actually observed will have a smaller range, as shown in Figure 9.1.

9.4 Mixture Melting Points

The melting point can be used as supporting evidence in identifying a compound in two different ways. Not only may the melting points of the two individual compounds be compared but a special procedure called a **mixture melting point** may also be performed. The mixture melting point requires that an authentic sample of the same compound be available from another source. In this procedure, the two compounds (authentic and suspected) are finely pulverized and mixed together in equal quantities. Then the melting point of the mixture is determined. If there is a melting-point depression or if the range of melting is expanded by a large amount compared to that of the individual substances, you may conclude that one compound has acted as an impurity toward the other and that they are not the same compound. If there is no lowering of the melting point for the mixture (the melting point is identical with those of pure A and pure B), then A and B are almost certainly the same compound.

9.5 Packing the Melting-Point Tube

Melting points are usually determined by heating the sample in a piece of thin-walled capillary tubing (1 mm × 100 mm) that has been sealed at one end. To pack the tube, press the open end gently into a *pulverized* sample of the crystalline material. Crystals will stick in the open end of the tube. The amount of solid pressed into the tube should correspond to a column no more than 1–2 mm high. To transfer the crystals to the closed end of the tube, drop the capillary tube, closed end first, down a $\frac{2}{3}$ -m length of glass tubing, which is held upright on the desktop. When the capillary tube hits the desktop, the crystals will pack down into the bottom of the tube. This procedure is repeated if necessary. Tapping the capillary on the desktop with fingers is not recommended because it is easy to drive the small tubing into a finger if the tubing should break.

Some commercial melting-point instruments have a built-in vibrating device that is designed to pack capillary tubes. With these instruments, the sample is pressed into the open end of the capillary tube, and the tube is placed in the vibrator slot. The action of the vibrator will transfer the sample to the bottom of the tube and pack it tightly.

9.6 Determining the Melting Point—the Thiele Tube

There are two principal types of melting-point apparatus available: the Thiele tube and commercially available, electrically heated instruments. The Thiele tube, shown in Figure 9.3, is the simpler device and was once widely used. It is a glass tube designed to contain a heating oil (mineral oil or silicone oil) and a thermometer to which a capillary tube containing the sample is attached. For best results

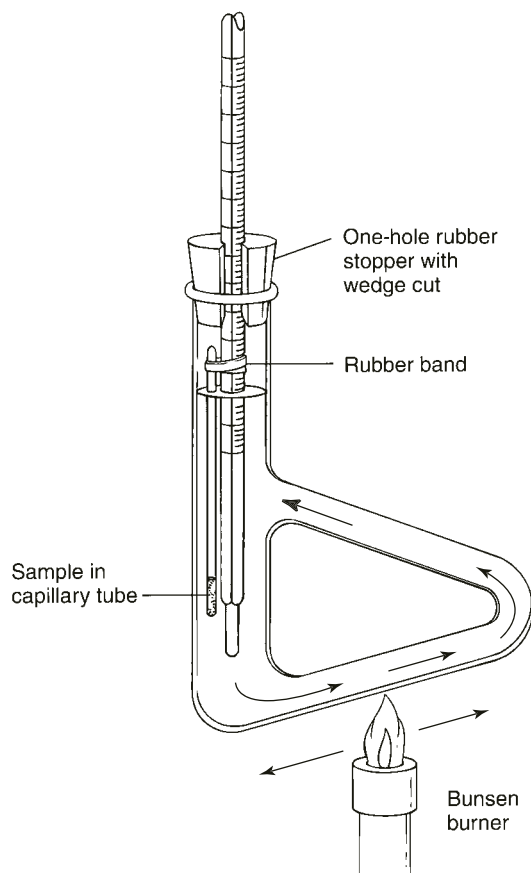


Figure 9.3
A Thiele tube.

use a mercury thermometer. The shape of the Thiele tube allows convection currents to form in the oil when it is heated. These currents maintain a uniform temperature distribution through the oil in the tube. The side arm of the tube is designed to generate these convection currents and thus transfer the heat from the flame evenly and rapidly throughout the oil. The sample, which is in a capillary tube attached to the thermometer, is held by a rubber band or a thin slice of rubber tubing. It is important that this rubber band be above the level of the oil (allowing for expansion of the oil on heating) so that the oil does not soften the rubber and allow the capillary tubing to fall into the oil. If a cork or a rubber stopper is used to hold the thermometer, a triangular channel should be cut into the side of it to allow pressure equalization.

The Thiele tube is usually heated by a microburner. During the heating, the rate of temperature increase should be regulated. Hold the burner by its cool base and, using a low flame, move the burner slowly back and forth along the bottom of the arm of the Thiele tube. If the heating is too fast, remove the burner for a few seconds and then resume heating. The rate of heating should be *slow* near the melting point (about 1°C per minute) to ensure that the temperature increase is not faster than the rate at which heat can be transferred to the sample being observed. At the melting point, it is necessary that the mercury in the thermometer and the sample in the capillary tube be at temperature equilibrium.

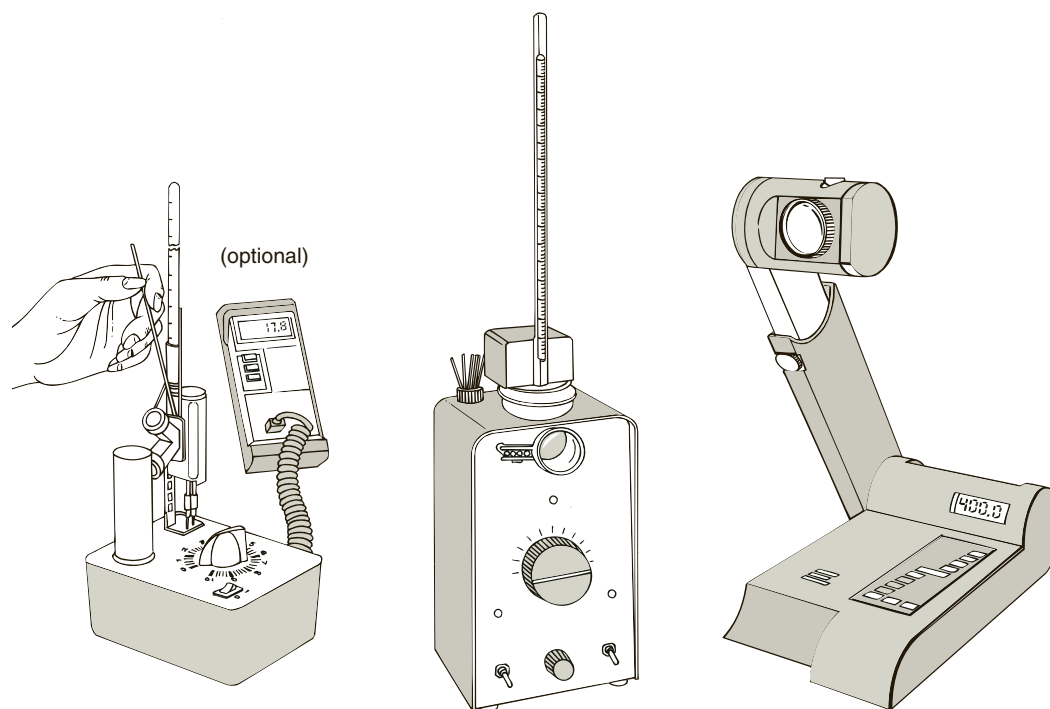


Figure 9.4
Melting-point apparatus.

9.7 Determining the Melting Point— Electrical Instruments

Three types of electrically heated melting-point instruments are illustrated in Figure 9.4. In each case, the melting-point tube is filled as described in Section 9.5 and placed in a holder located just behind the magnifying eyepiece. The apparatus is operated by moving the switch to the ON position, adjusting the potentiometric control dial for the desired rate of heating, and observing the sample through the magnifying eyepiece. The temperature is read from a thermometer or, in the most modern instruments, from a digital display attached to a thermocouple. Another option is the Vernier melt station, which can be connected to the Vernier Lab Quest unit or the Vernier LabPro system. (See Technique 13, Section 13.4 for more information about these devices.)

Your instructor will demonstrate and explain the type used in your laboratory.

Some electrically heated instruments do not heat or increase the temperature of the sample linearly. Although the rate of increase may be linear in the early stages of heating, it usually decreases and leads to a constant temperature at some upper limit. The upper-limit temperature is determined by the setting of the heating control. Thus, a family of heating curves is usually obtained for various control settings, as shown in Figure 9.5. The four hypothetical curves shown (1–4) might correspond to different control settings. For a compound melting at temperature t_1 , the setting corresponding to curve 3 would be ideal. In the beginning of the curve, the temperature is increasing too rapidly to allow determination of an

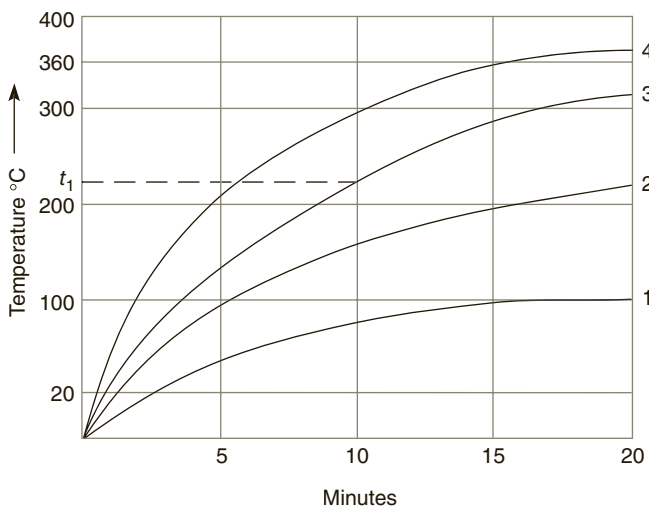


Figure 9.5
Heating-rate curves.

accurate melting point, but after the change in slope, the temperature increase will have slowed to a more usable rate.

If the melting point of the sample is unknown, you can often save time by preparing two samples for melting-point determination. With one sample, you can rapidly determine a crude melting-point value. Then repeat the experiment more carefully using the second sample. For the second determination, you already have an approximate idea of what the melting-point temperature should be, and a proper rate of heating can be chosen.

When temperatures above 150°C are measured, thermometer errors can become significant. For an accurate melting point with a high-melting solid, you may wish to apply a **stem correction** to the thermometer as described in Technique 13, Section 13.4. An even better solution is to calibrate the thermometer as described in Section 9.9.

9.8 Decomposition, Discoloration, Softening, Shrinkage, and Sublimation

Many solid substances undergo some degree of unusual behavior before melting. At times it may be difficult to distinguish these types of behavior from actual melting. You should learn, through experience, how to recognize melting and how to distinguish it from decomposition, discoloration, and particularly, softening and shrinkage.

Some compounds decompose on melting. This decomposition is usually evidenced by discoloration of the sample. Frequently, this decomposition point is a reliable physical property to be used in lieu of an actual melting point. Such decomposition points are indicated in tables of melting points by placing the symbol *d* immediately after the listed temperature. An example of a decomposition point is thiamine hydrochloride, whose melting point would be listed as 248°*d*, indicating that this substance melts with decomposition at 248°C. When decomposition is a result of reaction with the oxygen in air, it may be avoided by determining the melting point in a sealed, evacuated melting-point tube.

Figure 9.6 shows two simple methods of evacuating a packed tube. Method A uses an ordinary melting-point tube, and method B constructs the melting-point

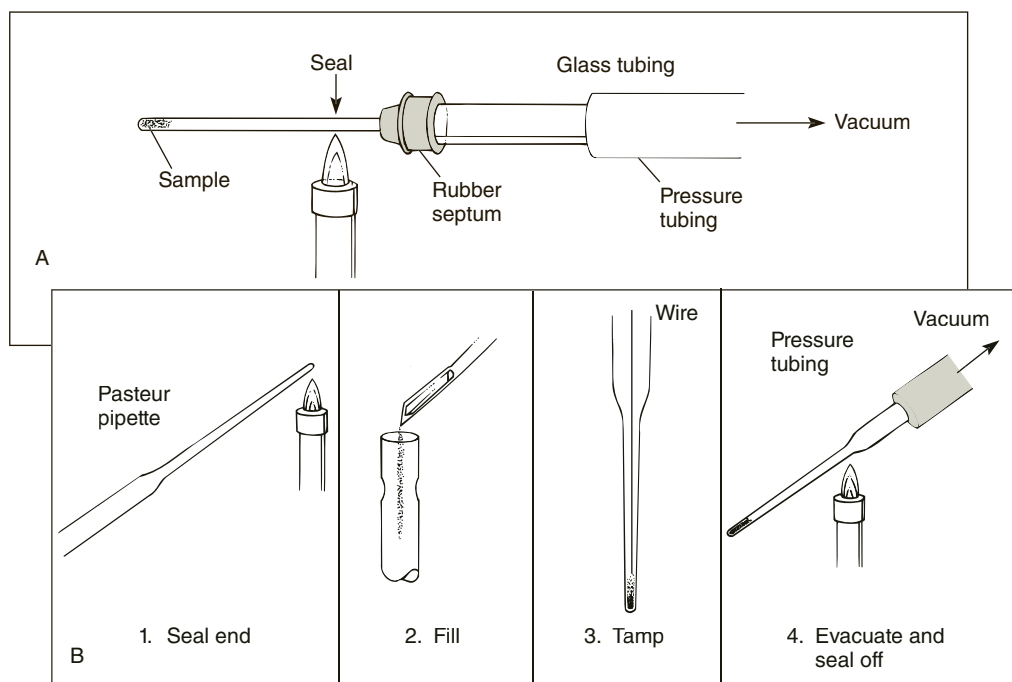


Figure 9.6
Evacuation and sealing of a melting-point capillary.

tube from a disposable Pasteur pipette. Before using method B, be sure to determine that the tip of the pipette will fit into the sample holder in your melting-point instrument.

Method A

In method A, a hole is punched through a rubber septum using a large pin or a small nail, and the capillary tube is inserted from the inside, sealed end first. The septum is placed over a piece of glass tubing connected to a vacuum line. After the tube is evacuated, the upper end of the tube may be sealed by heating and pulling it closed.

Method B

In method B, the thin section of a 9-inch Pasteur pipette is used to construct the melting-point tube. Carefully seal the tip of the pipette using a flame. Be sure to hold the tip *upward* as you seal it. This will prevent water vapor from condensing inside the pipette. When the sealed pipette has cooled, the sample may be added through the open end using a microspatula. A small wire may be used to compress the sample into the closed tip. (If your melting-point apparatus has a vibrator, it may be used in place of the wire to simplify the packing.) When the sample is in place, the pipette is connected to the vacuum line with tubing and evacuated. The evacuated sample tube is sealed by heating it with a flame and pulling it closed.

Some substances begin to decompose *below* their melting points. Thermally unstable substances may undergo elimination reactions or anhydride formation reactions during heating. The decomposition products formed represent impurities in the original sample, so the melting point of the substance may be lowered due to their presence.

It is normal for many compounds to soften or shrink immediately before melting. Such behavior represents not decomposition but a change in the crystal structure or a mixing with impurities. Some substances “sweat,” or release solvent of crystallization, before melting. These changes do not indicate the beginning of melting. Actual melting begins when the first drop of liquid becomes visible, and the melting range continues until the temperature is reached at which all the solid has been converted to the liquid state. With experience, you soon learn to distinguish between softening, or “sweating,” and actual melting. If you wish, the temperature of the onset of softening or sweating may be reported as a part of your melting-point range: 211°C (softens), 223–225°C (melts).

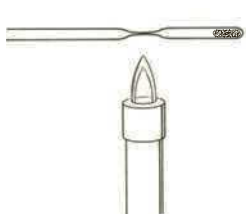


Figure 9.7
Sealing a tube for a substance that sublimes.

Some solid substances have such a high vapor pressure that they sublime at or below their melting points. In many handbooks, the sublimation temperature is listed along with the melting point. The symbols *sub*, *subl*, and sometimes *s* are used to designate a substance that sublimes. In such cases, the melting-point determination must be performed in a sealed capillary tube to avoid loss of the sample. The simplest way to seal a packed tube is to heat the open end of the tube in a flame and pull it closed with tweezers or forceps. A better way, although more difficult to master, is to heat the center of the tube in a small flame, rotating it about its axis, and keeping the tube straight, until the center collapses. If this is not done quickly, the sample may melt or sublime while you are working. With the smaller chamber, the sample will not be able to migrate to the cool top of the tube that may be above the viewing area. Figure 9.7 illustrates the method.

9.9 Thermometer Calibration

When a melting-point or boiling-point determination has been completed, you expect to obtain a result that duplicates the result recorded in a handbook or in

the original literature. It is not unusual, however, to find a discrepancy of several degrees from the literature value. Such a discrepancy does not necessarily indicate that the experiment was incorrectly performed or that the material is impure; rather, it may indicate that the thermometer used for the determination was slightly in error. Most thermometers do not measure the temperature with perfect accuracy.

To determine accurate values, you must calibrate the thermometer that is used. This calibration is done by determining the melting points of a variety of standard substances with the thermometer. A plot is drawn of the observed temperature vs. the published value of each standard substance. A smooth line is drawn through the points to complete the chart. A correction chart prepared in this way is shown in Figure 9.8. This chart is used to correct any melting point determined with that particular thermometer. Each thermometer requires its own calibration curve. A list of suitable standard substances for calibrating thermometers is provided in Table 9.1. The standard substances, of course, must be pure in order for the corrections to be valid.

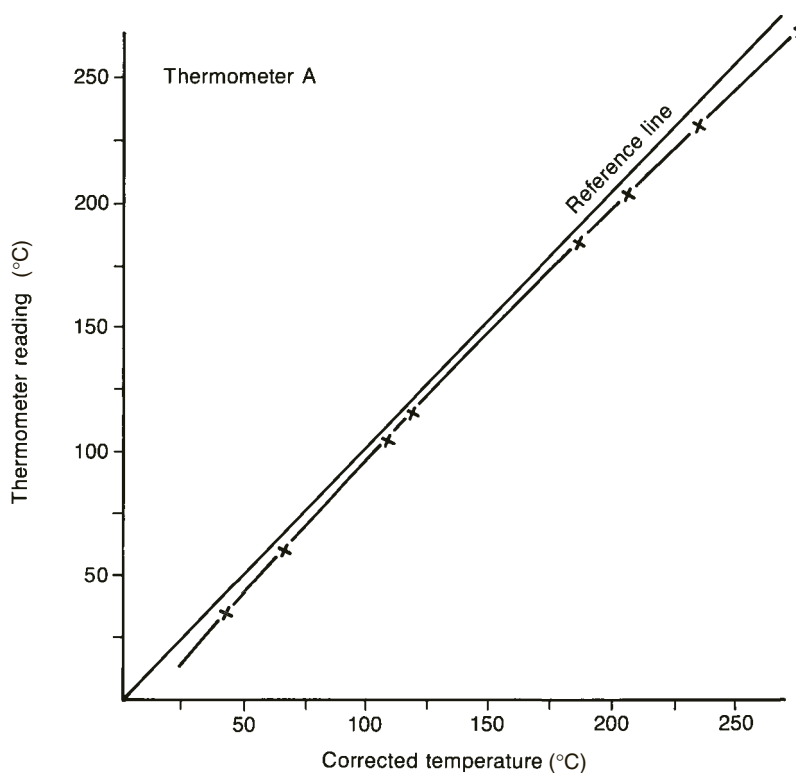


Figure 9.8
A thermometer-calibration curve.

Table 9.1 Melting-point standards

<i>Compound</i>	<i>Melting Point (°C)</i>
Ice (solid–liquid water)	0
Acetanilide	115
Benzamide	128
Urea	132
Succinic acid	189
3,5-Dinitrobenzoic acid	205

PROBLEMS

1. Two substances, A and B, have the same melting point. How can you determine if they are the same without using any form of spectroscopy? Explain in detail.
2. Using Figure 9.5, determine which heating curve would be most appropriate for a substance with a melting point of about 150°C.
3. What steps can you take to determine the melting point of a substance that sublimes before it melts?
4. A compound melting at 134°C was suspected to be either aspirin (mp 135°C) or urea (mp 133°C). Explain how you could determine whether one of these two suspected compounds was identical to the unknown compound without using any form of spectroscopy.
5. An unknown compound gave a melting point of 230°C. When the molten liquid solidified, the melting point was redetermined and found to be 131°C. Give a possible explanation for this discrepancy.

Solubility

The solubility of a **solute** (a dissolved substance) in a **solvent** (the dissolving medium) is the most important chemical principle underlying three basic techniques you will study in the organic chemistry laboratory: crystallization, extraction, and chromatography. In this discussion of solubility, you will gain an understanding of the structural features of a substance that determine its solubility in various solvents. This understanding will help you to predict solubility behavior and to understand the techniques that are based on this property. Understanding solubility behavior will also help you understand what is going on during a reaction, especially when there is more than one liquid phase present or when a precipitate is formed.

10.1 Definition of Solubility

Although we often describe solubility behavior in terms of a substance being **soluble** (dissolved) or **insoluble** (not dissolved) in a solvent, solubility can be described more precisely in terms of the *extent* to which a substance is soluble. Solubility may be expressed in terms of grams of solute per liter (g/L) or milligrams of solute per milliliter (mg/mL) of solvent. Consider the solubilities at room temperature for the following three substances in water:

Cholesterol	0.002 mg/mL
Caffeine	22 mg/mL
Citric acid	620 mg/mL

In a typical test for solubility, 40 mg of solute is added to 1 mL of solvent. Therefore, if you were testing the solubility of these three substances, cholesterol would be insoluble, caffeine would be partially soluble, and citric acid would be soluble. Note that a small amount (0.002 mg) of cholesterol would dissolve. It is unlikely, however, that you would be able to observe this small amount dissolving, and you would report that cholesterol is insoluble. On the other hand, 22 mg (55%) of the caffeine would dissolve. It is likely that you would be able to observe this, and you would state that caffeine is partially soluble.

When the solubility of a liquid solute in a solvent is described, it is sometimes helpful to use the terms **miscible** and **immiscible**. Two liquids that are miscible will mix homogeneously (one phase) in all proportions. For example, water and ethyl alcohol are miscible. When they are mixed in any proportion, only one layer will be observed. When two liquids are miscible, it is also true that either one of them will be completely soluble in the other one. Two immiscible liquids do not mix homogeneously in all proportions, and under some conditions they will form two layers. Water and diethyl ether are immiscible. When mixed in roughly equal amounts, they will form two layers. However, each liquid is slightly soluble in the other one. Even when two layers are present, a small amount of water will be soluble in the diethyl ether, and a small amount of diethyl ether will be soluble in the water. Furthermore, if only a small amount of either one is added to the other, it

may dissolve completely, and only one layer will be observed. For example, if a small amount of water (less than 1.2% at 20°C) is added to diethyl ether, the water will dissolve completely in the diethyl ether, and only one layer will be observed. When more water is added (more than 1.2%), some of the water will not dissolve, and two layers will be present.

Although the terms *solubility* and *miscibility* are related in meaning, it is important to understand that there is one essential difference. There can be different degrees of solubility, such as slightly, partially, very, and so on. Unlike solubility, miscibility does not have any degrees—a pair of liquids is either miscible or it is not.

10.2 Predicting Solubility Behavior

A major goal of this section is to explain how to predict whether a substance will be soluble in a given solvent. This is not always easy, even for an experienced chemist. However, guidelines will help you make a good guess about the solubility of a compound in a specific solvent. In discussing these guidelines, it is helpful to separate the types of solutions we will be looking at into two categories: solutions in which both the solvent and the solute are covalent (molecular), and ionic solutions, in which the solute ionizes and dissociates.

A. Solutions in Which the Solvent and Solute Are Molecular

A useful generalization in predicting solubility is the widely used rule “Like dissolves like.” This rule is most commonly applied to polar and nonpolar compounds. According to this rule, a polar solvent will dissolve polar (or ionic) compounds, and a nonpolar solvent will dissolve nonpolar compounds.

The reason for this behavior involves the nature of intermolecular forces of attraction. Although we will not be focusing on the nature of these forces, it is helpful to know what they are called. The force of attraction between polar molecules is called **dipole–dipole interaction**; between nonpolar molecules, forces of attraction are called **van der Waals forces** (also called **London** or **dispersion forces**). In both cases, these attractive forces can occur between molecules of the same compound or different compounds. Consult your lecture textbook for more information on these forces.

To apply the rule “Like dissolves like,” you must first determine whether a substance is polar or nonpolar. The polarity of a compound is dependent on both the polarities of the individual bonds and the shape of the molecule. For most organic compounds, evaluating these factors can become quite complicated because of the complexities of the molecules. However, it is possible to make some reasonable predictions just by looking at the types of atoms that a compound possesses. As you read the following guidelines, it is important to understand that although we often describe compounds as being polar or nonpolar, polarity is a matter of degree, ranging from nonpolar to highly polar.


Guidelines for Predicting Polarity and Solubility

1. All hydrocarbons are nonpolar.

Examples:



Hexane

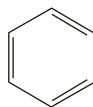


Benzene

Hydrocarbons such as benzene are slightly more polar than hexane because of their pi (π) bonds, which allow for greater van der Waals or London attractive forces.

2. Compounds possessing the electronegative elements oxygen or nitrogen are polar.


Examples:

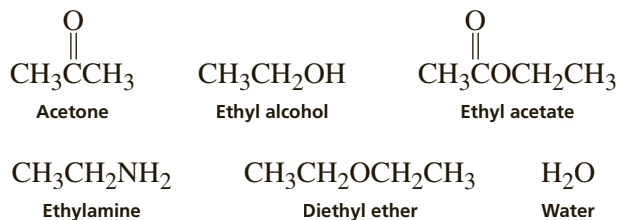


The polarity of these compounds depends on the presence of polar C—O, C=O, OH, NH, and CN bonds. The compounds that are most polar are capable of forming hydrogen bonds (see Guideline 6) and have NH or OH bonds. Although all these compounds are polar, the degree of polarity ranges from slightly polar to highly polar. This is due to the effect on polarity of the shape of the molecule and size of the carbon chain and whether the compound can form hydrogen bonds.

3. The presence of halogen atoms, even though their electronegativities are relatively high, does not alter the polarity of an organic compound in a significant way. Therefore, these compounds are only slightly polar. The polarities of these compounds are more similar to those of hydrocarbons, which are nonpolar, than to that of water, which is highly polar.

Examples:

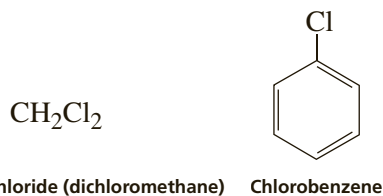

4. When comparing organic compounds within the same family, note that adding carbon atoms to the chain decreases the polarity. For example, methyl alcohol (CH_3OH) is more polar than propyl alcohol ($\text{CH}_3\text{CH}_2\text{CH}_2\text{OH}$). The reason is that hydrocarbons are nonpolar, and increasing the length of a carbon chain makes the compound more hydrocarbon-like.
5. Compounds that contain four or fewer carbons and also contain oxygen or nitrogen are often soluble in water. Almost any functional group containing these elements will lead to water solubility for low-molecular-weight (up to C_4) compounds. Compounds having five or six carbons and containing one of these elements are often insoluble in water or have borderline solubility.
6. As mentioned earlier, the force of attraction between polar molecules is dipole–dipole interaction. A special case of dipole–dipole interaction is hydrogen bonding. Hydrogen bonding is a possibility when a compound possesses a hydrogen atom bonded to a nitrogen, oxygen, or fluorine atom. The bond is formed by the attraction between this hydrogen atom and a nitrogen, oxygen, or fluorine atom in another molecule. Hydrogen bonding may occur between two molecules of the same compound or between molecules of different compounds:

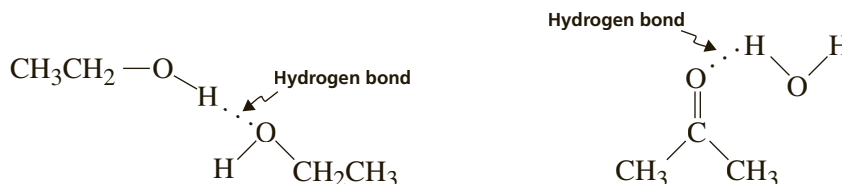

Hydrogen bonding is the strongest type of dipole–dipole interaction. When hydrogen bonding between solute and solvent is possible, solubility is greater than one would expect for compounds of similar polarity that cannot form hydrogen bonds. Hydrogen bonding is very important in organic chemistry, and you should be alert for situations in which hydrogen bonding may occur.

7. Another factor that can affect solubility is the degree of branching of the alkyl chain in a compound. Branching of the alkyl chain in a compound lowers the intermolecular forces between the molecules. This is usually reflected in a greater solubility in water for the branched compound than for the corresponding straight-chain compound. This occurs simply because the molecules of the branched compounds are more easily separated from one another.
8. The solubility rule (“Like dissolves like”) may be applied to organic compounds that belong to the same family. For example, 1-octanol (an alcohol) is soluble in the solvent ethyl alcohol. Most compounds within the same family have similar polarity. However, this generalization may not apply if there is a substantial difference in size between the two compounds. For example, cholesterol, an alcohol with a molecular weight (MW) of 386.64, is only slightly soluble in methanol (MW 32.04). The large hydrocarbon component of cholesterol negates the fact that they belong to the same family.
9. Almost all organic compounds that are in the ionic form are water soluble (see next section, B. Solutions in Which the Solute Ionizes and Dissociates).
10. The stability of the crystal lattice also affects solubility. Other things being equal, the higher the melting point (the more stable the crystal) is, the less soluble the compound. For instance, *p*-nitrobenzoic acid (mp 242°C) is, by a factor of 10, less soluble in a fixed amount of ethanol than the *ortho* (mp 147°C) and *meta* (mp 141°C) isomers.

You can check your understanding of some of these guidelines by studying the list given in Table 10.1, which is given in order of increasing polarity. The structures of these compounds are given on the preceding pages.

Table 10.1 Compounds in increasing order of polarity

	<i>Increasing Polarity</i>
Aliphatic hydrocarbons	
Hexane (nonpolar)	
Aromatic hydrocarbons (π bonds)	
Benzene (nonpolar)	
Halocarbons	
Methylene chloride (slightly polar)	
Compounds with polar bonds	
Diethyl ether (slightly polar)	
Ethyl acetate (intermediate polarity)	
Acetone (intermediate polarity)	
Compounds with polar bonds and hydrogen bonding	
Ethyl alcohol (intermediate polarity)	
Methyl alcohol (intermediate polarity)	
Water (highly polar)	

Table 10.2 Solvents in increasing order of polarity*Increasing Polarity (Approximate)*

RH	Alkanes (hexane, petroleum ether)
ArH	Aromatics (benzene, toluene)
ROR	Ethers (diethyl ether)
RX	Halides ($\text{CH}_2\text{Cl}_2 > \text{CHCl}_3 > \text{CCl}_4$)
RCOOR	Esters (ethyl acetate)
RCOR	Aldehydes, ketones (acetone)
RNH_2	Amines (triethylamine, pyridine)
ROH	Alcohols (methanol, ethanol)
RCONH_2	Amides (N, N-dimethylformamide)
RCOOH	Organic acids (acetic acid)
▼ H_2O	Water


This list can be used to make some predictions about solubility, based on the rule “Like dissolves like.” Substances that are close to one another on this list will have similar polarities. Thus, you would expect hexane to be soluble in methylene chloride but not in water. Acetone should be soluble in ethyl alcohol. On the other hand, you might predict that ethyl alcohol would be insoluble in hexane. However, ethyl alcohol is soluble in hexane because ethyl alcohol is somewhat less polar than methyl alcohol or water. This last example demonstrates that you must be careful in using the guidelines on polarity for predicting solubilities. Ultimately, solubility tests must be done to confirm predictions until you gain more experience.

The trend in polarities shown in Table 10.1 can be expanded by including more organic families. The list in Table 10.2 gives an approximate order for the increasing polarity of organic functional groups. It may appear that there are some discrepancies between the information provided in these two tables. The reason is that Table 10.1 provides information about specific compounds, whereas the trend shown in Table 10.2 is for major organic families and is approximate.

B. Solutions in Which the Solute Ionizes and Dissociates


Many ionic compounds are highly soluble in water because of the strong attraction between ions and the highly polar water molecules. This also applies to organic compounds that can exist as ions. For example, sodium acetate consists of Na^+ and CH_3COO^- ions, which are highly soluble in water. Although there are some exceptions, **you may assume that all organic compounds that are in the ionic form will be water soluble.**

The most common way by which organic compounds become ions is in acid–base reactions. For example, carboxylic acids can be converted to water-soluble salts when they react with dilute aqueous NaOH:

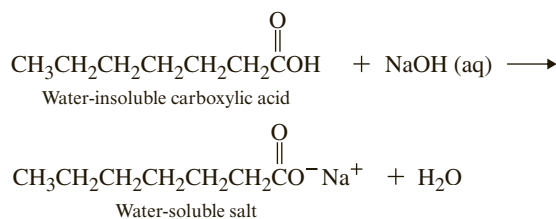

The water-soluble salt can then be converted back to the original carboxylic acid (which is insoluble in water) by adding another acid (usually aqueous HCl) to the solution of the salt. The carboxylic acid precipitates out of solution.

Amines, which are organic bases, can also be converted to water-soluble salts when they react with dilute aqueous HCl:

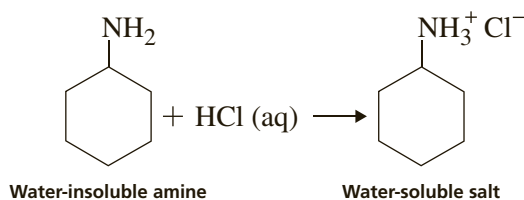

This salt can be converted back to the original amine by adding a base (usually aqueous NaOH) to the solution of the salt.

10.3 Organic Solvents

Organic solvents must be handled safely. Always remember that organic solvents are all at least mildly toxic and that many are flammable. You should become thoroughly familiar with laboratory safety (see Technique 1).

The most common organic solvents are listed in Table 10.3 along with their boiling points. Solvents marked in boldface type will burn. Ether, pentane, and hexane are especially dangerous; if they are combined with the correct amount of air, they will explode.

The terms **petroleum ether** and **ligroin** are often confusing. Petroleum ether is a mixture of hydrocarbons with isomers of formulas C_5H_{12} and C_6H_{14} predominating. Petroleum ether is not an ether at all because there are no oxygen-bearing compounds in the mixture. In organic chemistry, an ether is usually a compound containing an oxygen atom to which two alkyl groups are attached. Figure 10.1 shows some of the hydrocarbons that appear commonly in petroleum ether. It also

Table 10.3 Common organic solvents

<i>Solvent</i>	<i>Bp</i> (°C)	<i>Solvent</i>	<i>Bp</i> (°C)
Hydrocarbons		Ethers	
Pentane	36	Ether (diethyl)	35
Hexane	69	Dioxane ^a	101
Benzene ^a	80	1,2-Dimethoxyethane	83
Toluene	111	Others	
Hydrocarbon mixtures		Acetic acid	118
Petroleum ether	30–60	Acetic anhydride	140
Ligroin	60–90	Pyridine	115
Chlorocarbons		Acetone	56
Methylene chloride	40	Ethyl acetate	77
Chloroform ^a	61	Dimethylformamide	153
Carbon tetrachloride ^a	77	Dimethylsulfoxide	189
Alcohols			
Methanol	65		
Ethanol	78		
Isopropyl alcohol	82		

Note: **Boldface type** indicates flammability.

^aSuspected carcinogen.

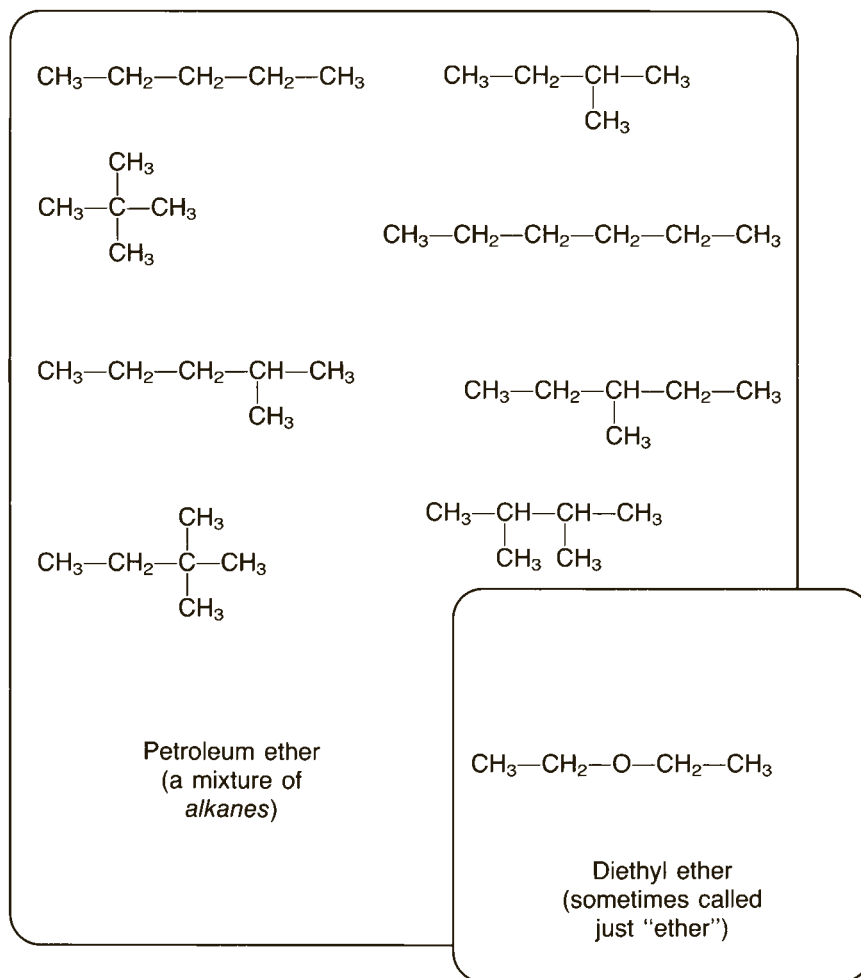


Figure 10.1
A comparison between “ether” (diethyl ether) and “petroleum ether.”

shows the structure of ether (diethyl ether). Use special care when instructions call for either **ether** or **petroleum ether**; the two must not become accidentally confused. Confusion is particularly easy when one is selecting a container of solvent from the supply shelf.


Ligroin, or high-boiling petroleum ether, is like petroleum ether in composition except that compared with petroleum ether, ligroin generally includes higher-boiling alkane isomers. Depending on the supplier, ligroin may have different boiling ranges. Whereas some brands of ligroin have boiling points ranging from about 60°C to about 90°C, other brands have boiling points ranging from about 60°C to about 75°C. The boiling-point ranges of petroleum ether and ligroin are often included on the labels of the containers.

PROBLEMS


- For each of the following pairs of solute and solvent, predict whether the solute would be soluble or insoluble. After making your predictions, you can check your answers by looking up the compounds in *The Merck Index* or the *CRC Handbook of Chemistry and Physics*. Generally, *The Merck Index* is the easier

reference book to use. If the substance has a solubility greater than 40 mg/mL, you may conclude that it is soluble.


a. Malic acid in water

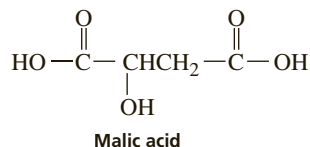


b. Naphthalene in water

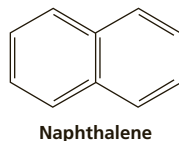


c. Amphetamine in ethyl alcohol

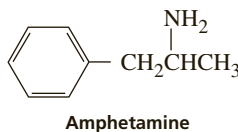


d. Aspirin in water

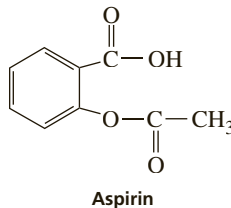


e. Succinic acid in hexane (*Note:* the polarity of hexane is similar to that of petroleum ether.)

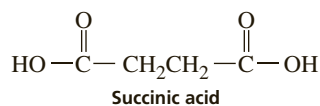

f. Ibuprofen in diethyl ether

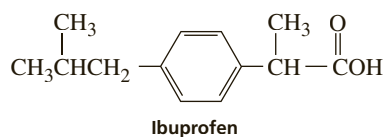

g. 1-Decanol (*n*-decyl alcohol) in water

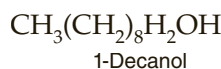


2. Predict whether the following pairs of liquids would be miscible or immiscible:

a. Water and methyl alcohol


b. Hexane and benzene

c. Methylene chloride and benzene


d. Water and toluene

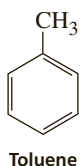


e. Ethyl alcohol and isopropyl alcohol

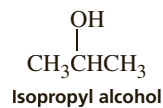

3. Would you expect ibuprofen (see problem 1f) to be soluble or insoluble in 1.0 M NaOH? Explain.
4. Thymol is very slightly soluble in water and very soluble in 1.0 M NaOH. Explain.

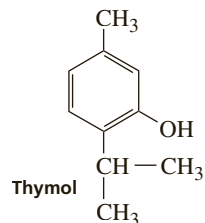

5. Although cannabinol and methyl alcohol are both alcohols, cannabinol is very slightly soluble in methyl alcohol at room temperature. Explain.

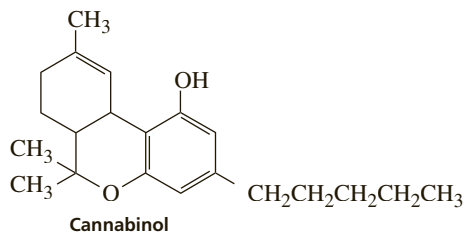

6. What is the difference between the compounds in each of the following pairs?
 - a. Ether and petroleum ether
 - b. Ether and diethyl ether
 - c. Ligroin and petroleum ether

Crystallization: Purification of Solids

In most organic chemistry experiments, the desired product is first isolated in an impure form. If this product is a solid, the most common method of purification is crystallization. The general technique involves dissolving the material to be crystallized in a *hot* solvent (or solvent mixture) and cooling the solution slowly. The dissolved material has a decreased solubility at lower temperatures and will separate from the solution as it is cooled. This phenomenon is called either **crystallization**, if the crystal growth is relatively slow and selective, or **precipitation**, if the process is rapid and nonselective. Crystallization is an equilibrium process and produces very pure material. A small seed crystal is formed initially, and it then grows layer by layer in a reversible manner. In a sense, the crystal “selects” the correct molecules from the solution. In precipitation, the crystal lattice is formed so rapidly that impurities are trapped within the lattice. Therefore, any attempt at purification with too rapid a process should be avoided. Because the impurities are usually present in much smaller amounts than the compound being crystallized, most of the impurities will remain in the solvent even when it is cooled. The purified substance can then be separated from the solvent and from the impurities by filtration.

In microscale organic work, two methods are commonly used to perform crystallizations. The first method, which is carried out with an Erlenmeyer flask to dissolve the material and a Hirsch funnel to filter the crystals, is normally used when the weight of solid to be crystallized is more than 0.1 g. This technique, called **semi-microscale crystallization**, is discussed in Section 11.3. The second method is performed with a Craig tube and is used with smaller amounts of solid. Referred to as **microscale crystallization**, this technique is discussed in Section 11.4. The weight of solid to be crystallized, however, is not the only factor to consider when choosing a method for crystallization. Because the solubility of a substance in a given solvent must also be taken into account, the weight, 0.1 g, should not be adhered to rigidly in determining which method to use. In this textbook, you will usually be advised which method to use in the experimental procedure.

The method described here for semimicroscale crystallizations is nearly identical to that used for crystallizing larger amounts of materials than those encountered in this textbook. Therefore, this technique can also be used to perform crystallizations at the macroscale level (more than several grams).

PART A. THEORY

11.1 Solubility

The first problem in performing a crystallization is selecting a solvent in which the material to be crystallized shows the desired solubility behavior. In an ideal case, the material should be sparingly soluble at room temperature and yet quite soluble

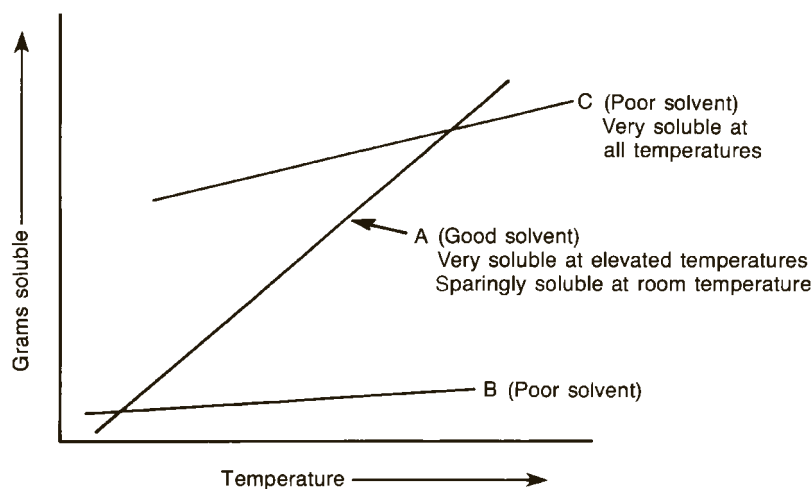


Figure 11.1
Graph of solubility vs. temperature.

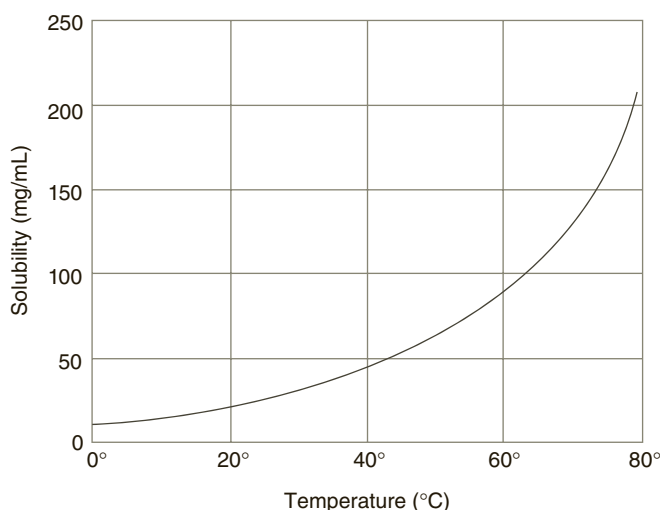


Figure 11.2
Solubility of sulfanilamide in 95% ethyl alcohol.

at the boiling point of the solvent selected. The solubility curve should be steep, as can be seen in line A of Figure 11.1. A curve with a low slope (line B) would not cause significant crystallization when the temperature of the solution was lowered. A solvent in which the material is very soluble at all temperatures (line C) also would not be a suitable crystallization solvent. The basic problem in performing a crystallization is to select a solvent (or mixed solvent) that provides a steep solubility-vs.-temperature curve for the material to be crystallized. A solvent that allows the behavior shown in line A is an ideal crystallization solvent. It should also be mentioned that solubility curves are not always linear, as they are depicted in Figure 11.1. This figure represents an idealized form of solubility behavior. The solubility curve for sulfanilamide in 95% ethyl alcohol, shown in

Figure 11.2, is typical of many organic compounds and shows what solubility behavior might look like for a real substance. This graph is based on the data in the following table:

<u>Temperature</u>	<u>Solubility (mg/mL)</u>
0°C	14
20°C	24
40°C	46
60°C	88
80°C	210

The solubility of organic compounds is a function of the polarities of both the solvent and the **solute** (dissolved material). A general rule is “Like dissolves like.” If the solute is very polar, a very polar solvent is needed to dissolve it; if the solute is nonpolar, a nonpolar solvent is needed. Applications of this rule are discussed extensively in Technique 10, Section 10.2, and in Section 11.5.

11.2 Theory of Crystallization

A successful crystallization depends on a large difference between the solubility of a material in a hot solvent and its solubility in the same solvent when it is cold. When the impurities in a substance are equally soluble in both the hot and the cold solvent, an effective purification is not easily achieved through crystallization. A material can be purified by crystallization when both the desired substance and the impurity have similar solubilities, but only when the impurity represents a small fraction of the total solid. The desired substance will crystallize on cooling, but the impurities will not.

For example, consider a case in which the solubilities of substance A and its impurity B are both 1 g/100 mL of solvent at 20°C and 10 g/100 mL of solvent at 100°C. In the impure sample of A, the composition is 9 g of A and 2 g of B. In the calculations for this example, it is assumed that the solubilities of both A and B are unaffected by the presence of the other substance. To make the calculations easier to understand, 100 mL of solvent are used in each crystallization. Normally, the minimum amount of solvent required to dissolve the solid would be used.

At 20°C, this total amount of material will not dissolve in 100 mL of solvent. However, if the solvent is heated to 100°C, all 11 g dissolve. The solvent has the capacity to dissolve 10 g of A and 10 g of B at this temperature. If the solution is cooled to 20°C, only 1 g of each solute can remain dissolved, so 8 g of A and 1 g of B crystallize, leaving 2 g of material in the solution. This crystallization is shown in Figure 11.3. The solution that remains after a crystallization is called the **mother liquor**. If the process is now repeated by treating the crystals with 100 mL of fresh solvent, 7 g of A will crystallize again, leaving 1 g of A and 1 g of B in the mother liquor. As a result of these operations, 7 g of pure A are obtained, but with the loss of 4 g of material (2 g of A plus 2 g of B). Again, this second crystallization step is illustrated in Figure 11.3. The final result illustrates an important aspect of crystallization—it is wasteful. Nothing can be done to prevent this waste; some A must be lost along with the impurity B for the method to be successful. Of course, if the impurity B were *more* soluble than A in the solvent, the losses would be reduced. Losses could also be reduced if the impurity were present in *much smaller* amounts than the desired material.

Note that in the preceding case, the method operated successfully because A was present in substantially larger quantity than its impurity B. If there had been a 50-50 mixture of A and B initially, no separation would have been achieved. In general, a crystallization is successful only if there is a *small* amount of impurity. As the amount of impurity increases, the loss of material must also increase. Two substances with nearly equal solubility behavior, present in equal amounts, cannot be separated. If the solubility behavior of two components present in equal amounts is different, however, a separation or purification is frequently possible.

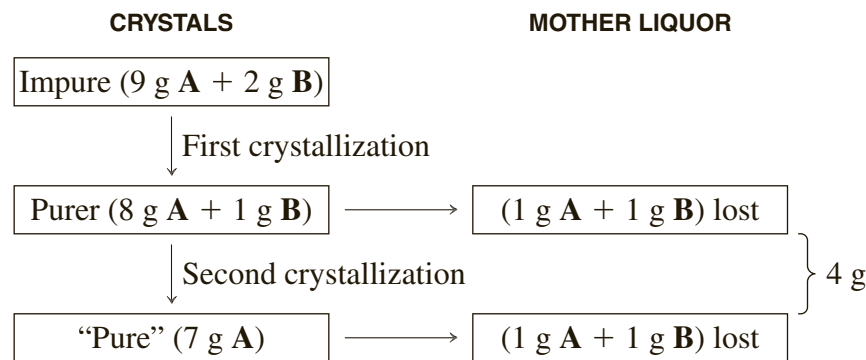


Figure 11.3
Purification of a mixture by crystallization.

In the preceding example, two crystallization procedures were performed. Normally, this is not necessary; however, when it is, the second crystallization is more appropriately called **recrystallization**. As illustrated in this example, a second crystallization results in purer crystals, but the yield is lower.

In most experiments, you will cool the crystallizing mixture in an ice-water bath before collecting the crystals by filtration. Cooling the mixture increases the yield by decreasing the solubility of the substance; however, even at this reduced temperature, some of the product will be soluble in the solvent. It is not possible to recover all your product in a crystallization procedure even when the mixture is cooled in an ice-water bath. A good example of this is illustrated by the solubility curve for sulfanilamide shown in Figure 11.2. The solubility of sulfanilamide at 0°C is still significant, 14 mg/mL. You should also remember that cooling the crystallizing mixture in an ice-water bath may cause more impurities to come out of solution, too, and the crystals may be less pure.

PART B. SEMIMICROSCALE CRYSTALLIZATION

11.3 Semimicroscale Crystallization—Hirsch Funnel

The crystallization technique described in this section is used when the weight of solid to be crystallized is more than 0.1 g. The four main steps in a semimicroscale crystallization are

1. Dissolving the solid
2. Removing insoluble impurities (when necessary)
3. Crystallization
4. Isolation of crystals

These steps are illustrated in Figure 11.4. It should be pointed out that a microscale crystallization with a Craig tube involves the same four steps, although the apparatus and procedures are somewhat different (see Section 11.4).

A. Dissolving the Solid

To minimize losses of material to the mother liquor, it is desirable to *saturate* the boiling solvent with solute. This solution, when cooled, will return the maximum possible amount of solute as crystals. To achieve this high return, the solvent is brought to its boiling point, and the solute is dissolved in the *minimum amount (!) of boiling solvent*. For this procedure, it is advisable to maintain a container of boiling solvent (on either a hot plate or a sand bath). From this container, add a small portion (about 0.5 mL) of the solvent to the flask¹ (usually a 10- or 25-mL Erlenmeyer flask) containing the solid to be crystallized and heat this mixture while swirling occasionally until it resumes boiling.


CAUTION



Do not heat the flask containing the solid until after you have added the first portion of solvent.

If the solid does not dissolve in the first portion of boiling solvent, then another small portion of boiling solvent is added to the flask. The mixture is heated again until

¹ A beaker should not be used, because the large opening allows the solvent to evaporate too rapidly and dust particles to get in too easily.

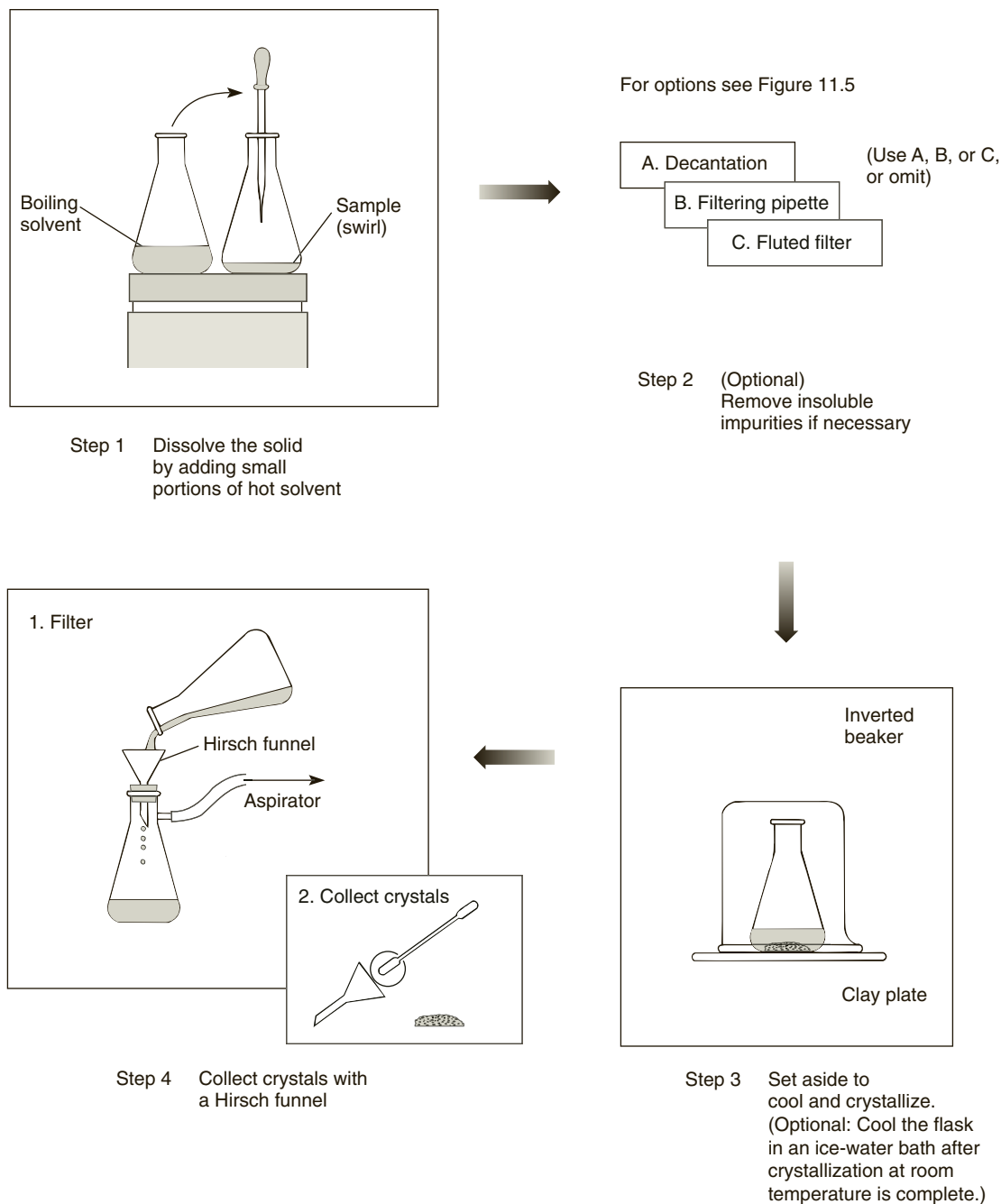**Figure 11.4**

Steps in a semimicroscale crystallization (no decolorization).

it resumes boiling. If the solid dissolves, no more solvent is added. But if the solid has not dissolved, another portion of boiling solvent is added as before, and the process is repeated until the solid dissolves. (If the solid totally dissolves in less than 2 mL of solvent, a Craig tube should be used for crystallization.) The portions of solvent added each time should be small so that only the *minimum* amount of solvent necessary for dissolving the solid is added. It is also important to emphasize that the procedure requires the addition of solvent to solid. You must never add portions of solid to a fixed quantity of boiling solvent. By this latter method, it is impossible to tell when saturation

has been achieved. You should perform this procedure fairly rapidly. Otherwise, you may lose solvent by evaporation nearly as quickly as you are adding it, and dissolving the solid will take a long time. This is most likely to happen when highly volatile solvents such as methyl alcohol or ethyl alcohol are used. The time from the first addition of solvent until the solid dissolves completely should not be longer than 10–15 minutes.

Comments on This Procedure for Dissolving the Solid

1. One of the most common mistakes is to add too much solvent. This can happen most easily if the solvent is not hot enough or if the mixture is not stirred sufficiently. If too much solvent is added, the percentage recovery will be reduced; it is even possible that no crystals will form when the solution is cooled. If too much solvent is added, you must evaporate the excess by heating the mixture. A nitrogen or air stream directed into the container will accelerate the evaporation process (see Technique 7, Section 7.10).
2. Another common mistake is dissolving the solid in hot solvent that is *below* the boiling point of the solvent. This can significantly reduce the recovery.
3. It is very important not to heat the solid until you have added some solvent. Otherwise, the solid may melt and possibly form an oil or decompose, and it may not crystallize easily (see Section 11.5).
4. It is also important to use an Erlenmeyer flask rather than a beaker for performing the crystallization. A beaker should not be used because the large opening allows the solvent to evaporate too rapidly and allows dust particles to get in too easily.
5. In some experiments, a specified amount of solvent for a given weight of solid will be recommended. In these cases, you should use the amount specified rather than the minimum amount of solvent necessary to dissolve the solid. The amount of solvent recommended has been selected to provide the optimum conditions for good crystal formation.
6. Occasionally, you may encounter an impure solid that contains small particles of insoluble impurities, pieces of dust, or paper fibers that will not dissolve in the hot crystallizing solvent. A common error is to add too much of the hot solvent in an attempt to dissolve these small particles, not realizing that they are insoluble. In such cases, you must be careful not to add too much solvent.
7. It is sometimes necessary to decolorize the solution by adding activated charcoal or by passing the solution through a column containing alumina or silica gel (see Section 11.7 and Technique 19, Section 19.14). A decolorization step should be performed only if the mixture is *highly* colored and it is clear that the color is due to impurities and not to the actual color of the substance being crystallized. If decolorization is necessary, it should be accomplished before the following filtration step.

B. Removing Insoluble Impurities

It is necessary to use one of the following three methods only if insoluble material remains in the hot solution or if decolorizing charcoal has been used.

NOTE: Indiscriminate use of the procedure can lead to needless loss of your product.

Decantation is the easiest method of removing solid impurities and should be considered first. A filtering pipette is used when the volume of liquid to be filtered is less than 10 mL (see Technique 8, Section 8.1, Part C), and you should use gravity

filtration through a fluted filter when the volume is 10 mL or greater (see Technique 8, Section 8.1, Part B). These three methods are illustrated in Figure 11.5.

Decantation. If the solid particles are relatively large or they easily settle to the bottom of the flask, it may be possible to separate the hot solution from the impurities by carefully pouring off the liquid, leaving the solid behind. This is accomplished most easily by holding a glass stirring rod along the top of the flask and tilting the flask so that the liquid pours out along one end of the glass rod into another container. A technique similar in principle to decantation, which may be easier to perform with smaller amounts of liquid, is to use a **preheated Pasteur pipette** to remove the hot solution. With this method, it may be helpful to place the tip of the pipette against the bottom of the flask when removing the last portion of solution. The small space between the tip of the pipette and the inside surface of the flask prevents solid material from being drawn into the pipette. An easy way to preheat the pipette is to draw up a small portion of hot *solvent* (not the *solution* being transferred) into the pipette and expel the liquid. Repeat this process several times.

Filtering Pipette. If the volume of solution after dissolving the solid in hot solvent is less than 10 mL, gravity filtration with a filtering pipette may be used to remove solid impurities. However, using a filtering pipette to filter a hot solution saturated with solute can be difficult without premature crystallization. The best way to prevent this from occurring is to add enough solvent to dissolve the desired product at room temperature (be sure not to add too much solvent) and perform the filtration at room temperature, as described in Technique 8, Section 8.1, Part C. After filtration, the excess solvent is evaporated by boiling until the solution is saturated at the boiling point of the mixture (see Technique 7, Section 7.10). If powdered decolorizing charcoal was used, it will probably be necessary to perform two filtrations with a filtering pipette, or else the method described next can be used.

Fluted Filter. This method is the most effective way to remove solid impurities when the volume of liquid is greater than 10 mL or when decolorizing charcoal has been used (see Technique 8, Section 8.1, Part B). You should add a small amount of extra solvent to the hot mixture. This procedure helps to prevent crystal formation in the filter paper or the stem of the funnel during the filtration. The funnel is fitted with a fluted filter and installed at the top of the Erlenmeyer flask to be used for the actual filtration. It is advisable to place a small piece of wire between the funnel and the mouth of the flask to relieve any increase in pressure caused by hot filtrate.

The Erlenmeyer flask containing the funnel and the fluted paper is placed on top of a sand bath or hot plate (low setting). The liquid to be filtered is brought to its boiling point and poured through the filter in portions. (If the volume of the mixture is less than 10 mL, it may be more convenient to transfer the mixture to the filter with a preheated Pasteur pipette.) It is necessary to keep the solutions in both flasks at their boiling temperatures to prevent premature crystallization. The refluxing action of the filtrate keeps the funnel warm and reduces the chance that the filter will clog with crystals that may have formed during the filtration. With low-boiling solvents, be aware that some solvent may be lost through evaporation. Consequently, extra solvent must be added to make up for this loss. If crystals begin to form in the filter during filtration, a minimum amount of boiling solvent is added to redissolve the crystals and to allow the solution to pass through the funnel. If the volume of liquid being filtered is less than 10 mL, a small amount of hot solvent should be used to rinse the filter after all the filtrate has been collected. The rinse solvent is then combined with the original filtrate.

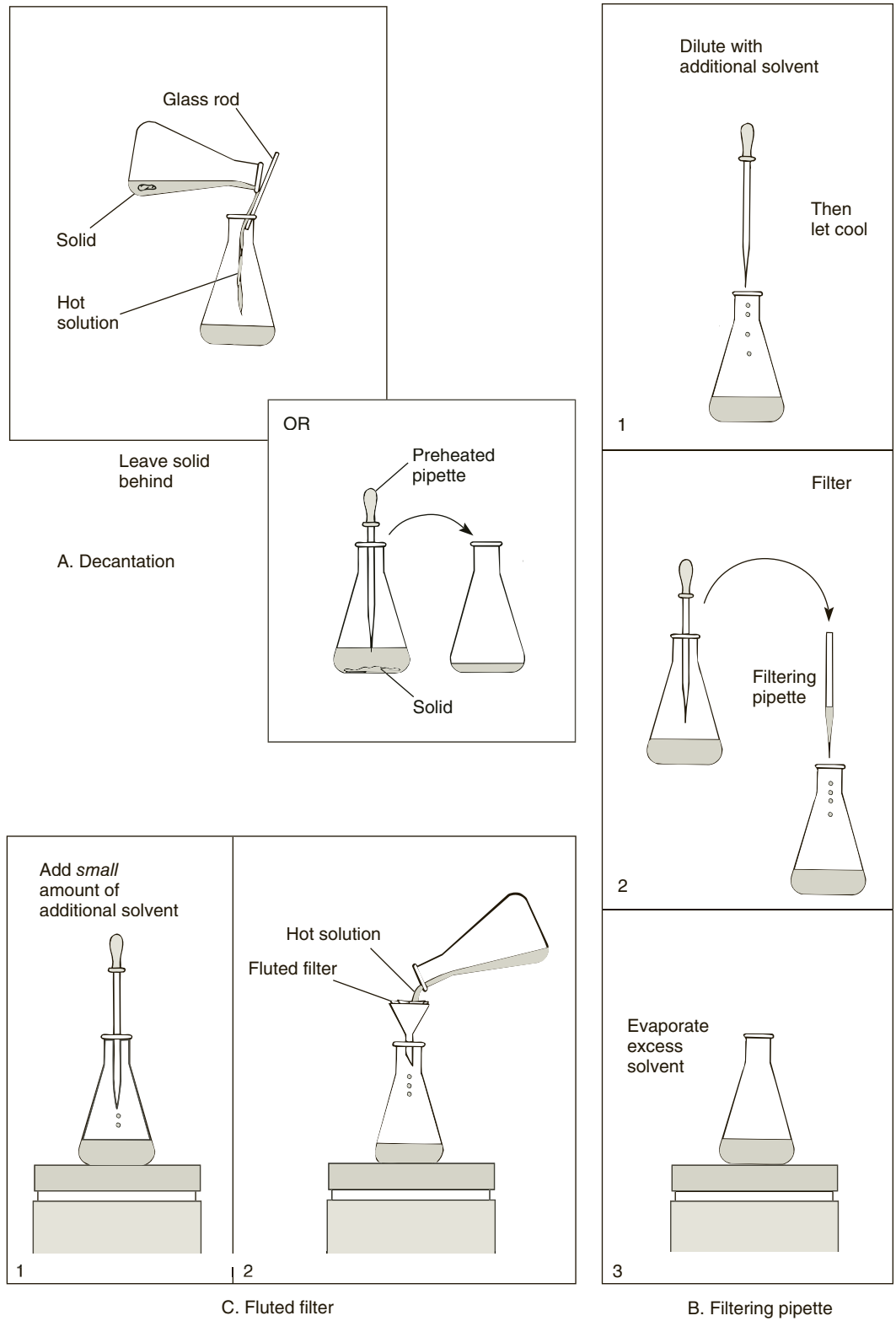


Figure 11.5
Methods for removing insoluble impurities in a semimicroscale crystallization.

After the filtration, it may be necessary to remove extra solvent by evaporation until the solution is saturated at the boiling point of the solvent (see Technique 7, Section 7.10).

C. Crystallizing the Solid

An Erlenmeyer flask, not a beaker, should be used for crystallization. The large open top of a beaker makes it an excellent dust catcher. The narrow opening of the Erlenmeyer flask reduces contamination by dust and allows the flask to be stoppered if it is to be set aside for a long period. Mixtures set aside for long periods must be stoppered after cooling to room temperature to prevent evaporation of solvent. If all the solvent evaporates, no purification is achieved, and the crystals originally formed become coated with the dried contents of the mother liquor. Even if the time required for crystallization to occur is relatively short, it is advisable to cover the top of the Erlenmeyer flask with a small watch glass or inverted beaker to prevent evaporation of solvent while the solution is cooling to room temperature.

The chances of obtaining pure crystals are improved if the solution cools to room temperature slowly. When the volume of solution is 10 mL or less, the solution is likely to cool more rapidly than is desired. This can be prevented by placing the flask on a surface that is a poor heat conductor and covering the flask with a beaker to provide a layer of insulating air. Appropriate surfaces include a clay plate or several pieces of filter paper on top of the laboratory bench. It may also be helpful to use a clay plate that has been warmed slightly on a hot plate or in an oven.

After crystallization at room temperature is complete, it is usually desirable to cool the flask in an ice-water bath. Because the solute is less soluble at lower temperatures, this will increase the yield of crystals.

If a cooled solution does not crystallize, it will be necessary to induce crystallization. Several techniques are described in Section 11.8, Part A.

D. Collecting and Drying the Crystals

After the flask has been cooled, the crystals are collected by vacuum filtration through a Hirsch (or Büchner) funnel (see Technique 8, Section 8.3, and Figure 8.5).

Moisten the filter paper with a few drops of the crystallizing solvent and turn on the vacuum (or aspirator) to the fullest extent. Use a spatula to dislodge the crystals from the bottom of the flask before transferring the material to the funnel. Swirl the mixture in the flask and pour the mixture into the funnel, attempting to transfer both crystals and solvent. You need to pour the mixture quickly, before the crystals have completely resettled on the bottom of the flask. You may need to do this in portions, depending on the size of your funnel. If all of the liquid and crystals will fit into the funnel, the most effective way to make this transfer is to “swirl and dump” all of the mixture as quickly as you can. By doing this, you are more likely to transfer all or most of the crystals without needing to add more solvent to transfer the remaining crystals.

The crystals should be washed with a small amount of *cold* solvent to remove any mother liquor adhering to their surface. Hot or warm solvent will dissolve some of the crystals. The crystals should then be left for a short time in the funnel, where air, as it passes, will dry them free of most of the solvent. It is often wise to cover the Hirsch funnel with an oversize filter paper or towel during this air-drying. This precaution prevents accumulation of dust in the crystals. When the crystals are nearly dry, they should be gently scraped off (so paper fibers are not removed with the crystals) the filter paper onto a watch glass or clay plate for further drying (see Section 11.9).

PART C. MICROSCALE CRYSTALLIZATION

11.4 Microscale Crystallization— Craig Tube

In most microscale experiments, the amount of solid to be crystallized is small enough (generally less than 0.1 g) that a **Craig tube** (see Technique 8, Figure 8.10) is the preferred method for crystallization. The main advantage of the Craig tube is that it minimizes the number of transfers of solid material, thus resulting in a greater yield of crystals. Also, the separation of the crystals from the mother liquor with the Craig tube is very efficient, and little time is required for drying the crystals. The steps involved are, in principle, the same as those performed when a crystallization is accomplished with an Erlenmeyer flask and a Hirsch funnel. The steps in a microscale crystallization using a Craig tube are illustrated in Figure 11.6.

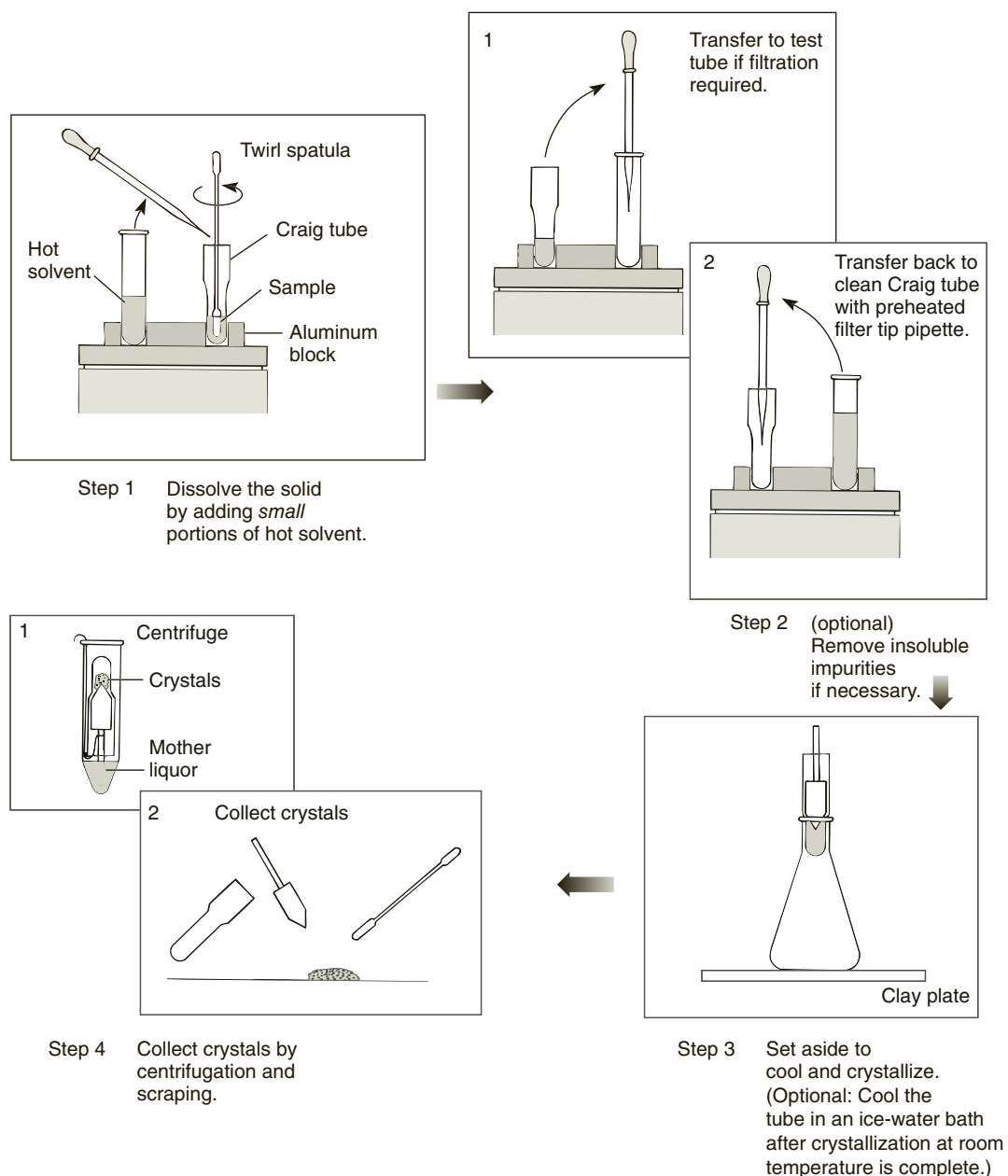


Figure 11.6
Steps in a microscale crystallization (no decolorization).

A. Dissolving the Solid

In crystallizations in which a filtration step is not required to remove insoluble impurities such as dirt or activated charcoal, this first step can be performed directly in the Craig tube. Otherwise, use a small test tube. The solid is transferred to the Craig tube, and the appropriate solvent, contained in a test tube, is heated to boiling on an aluminum block. A small portion (several drops) of hot solvent is added to the Craig tube, which is subsequently heated on the aluminum block until the solution in the Craig tube starts to boil. The hot mixture should be stirred continuously with a microspatula using a twirling motion. Stirring not only helps to dissolve the solute but also prevents the boiling liquid from bumping. Additional portions of hot solvent are added until all the solid has dissolved. In order to obtain the maximum yield, it is important not to add too much solvent, although any excess solvent can be evaporated later. You should perform this procedure fairly rapidly. Otherwise you may lose solvent by evaporation nearly as quickly as you are adding it, and dissolving all the solid will take a long time. The time required to dissolve the solid should not be longer than 15 minutes.

In many of the experiments in this textbook, a specified amount of solvent for a given weight of solid is recommended. In these cases, use the amount specified rather than the minimum amount of solvent necessary for dissolving the solid. The amount of solvent recommended has been selected to provide the optimum conditions for good crystal formation.

If the mixture is *highly* colored, and it is clear that the color is due to impurities and not to the actual color of the substance being crystallized, it will be necessary to decolorize the liquid. If decolorization is necessary, it should be accomplished before the following filtration step. Decolorizing charcoal may be used or the mixture may be passed through an alumina or silica gel column (see Section 11.7, Parts B and C, and Technique 19, Section 19.14).

B. Removing Insoluble Impurities

You should be alert for the presence of impurities that will not dissolve in the hot solvent, no matter how much solvent is added. If it appears that most of the solid has dissolved and the remaining solid has no tendency to dissolve, or if the liquid has been decolorized with charcoal, it will be necessary to remove the solid particles. Two methods are discussed.

If the impurities are relatively large or concentrated in one part of the mixture, it may be possible to use a Pasteur pipette preheated with hot solvent to draw up the liquid without removing any solid. One way to do this is to expel the air from the pipette and then place the end of the pipette on the bottom of the tube, being careful not to trap any solid in the pipette. The small space between the pipette and the bottom of the tube should allow you to draw up the liquid without removing any solid.

When filtration is necessary, a preheated Pasteur pipette is used to transfer the mixture to a test tube. After this transfer is made the Craig tube is rinsed with a few drops of solvent, which are also added to the test tube. The Craig tube is then washed and dried. The test tube containing the mixture is also heated in the sand bath. An additional 5 to 10 drops of solvent are added to the test tube to ensure that premature crystallization does not occur during the filtration step. To filter the mixture, take up the mixture in a filter-tip pipette (see Technique 8, Section 8.6) that has been preheated with hot solvent and quickly transfer the liquid to the clean Craig tube. Passing the liquid through the cotton plug in the filter-tip pipette should remove the

solid impurities. If this does not occur, it may be necessary to add more solvent (to prevent crystallization) and filter the mixture through a filtering pipette (Technique 8, Section 8.1, Part C). In either case, once the filtered solution has been returned to the Craig tube, it will be necessary to evaporate some solvent until the solution is saturated near the boiling point of the liquid. This is most conveniently accomplished by placing the Craig tube in the sand bath and, while stirring rapidly using a microspatula (twirling is most effective), bringing the solution to a boil. When you begin to observe a trace of solid material coating the spatula just above the level of the liquid, the solution is near saturation, and evaporation should be stopped.

C. Crystallizing the Solid

The hot solution is cooled slowly in the Craig tube to room temperature. Recall that slow cooling is important in the formation of pure crystals. When the volume of solution is 2 mL or less and the mass of glassware is relatively small, slow cooling is somewhat difficult to achieve. One method of increasing the cooling time is to insert the inner plug into the outer part of the Craig tube and place the Craig tube into a 10-mL Erlenmeyer flask. The layer of air in the flask will help insulate the hot solution as it cools. The Erlenmeyer flask is placed on a surface such as a clay plate (warmed slightly, if desired) or several pieces of paper. Another method is to fill a 10-mL Erlenmeyer flask with 8–10 mL of hot water at a temperature below the boiling point of the solvent. The assembled Craig tube is placed in the Erlenmeyer flask, which is set on an appropriate surface. Be careful not to put so much water in the Erlenmeyer flask that the Craig tube floats. After crystallization at room temperature is complete, the Craig tube can be placed in an ice-water bath to maximize the yield.

If crystals have not formed after the solution has cooled, it will be necessary to induce crystallization. Several techniques are described in Section 11.8.

A common occurrence with crystallizations using a Craig tube is to obtain a seemingly solid mass of small crystals. This may not be a problem, but if there is very little mother liquor present or the crystals are impure, it may be necessary to repeat the crystallization. This situation may have resulted either because the cooling process occurred too rapidly or because the solubility–temperature curve was so steep for a given solvent that very little mother liquor remained after the crystallization. In either case, you may want to repeat the crystallization to obtain a better (purer) yield of crystals. Three measures may be taken to avoid this problem. A small amount of extra solvent may be added before heating the mixture again and allowing it to cool. A second measure is to cool the solution more slowly. Finally, it may be helpful to try to induce crystallization *before* the solution has cooled to room temperature.

D. Collecting and Drying the Crystals

When the crystals have formed and the mixture has cooled in an ice-water bath (if desired), the Craig tube is placed in a centrifuge tube and the crystals are separated from the mother liquor by centrifugation (see Technique 8, Section 8.8). The crystals are then scraped off the end of the inner plug or from inside the Craig tube onto a watch glass or piece of paper. Minimal drying will be necessary (see Section 11.9).

The four steps in a semimicroscale or microscale crystallization are summarized in Table 11.1.

Table 11.1 Steps in a semimicroscale or microscale crystallization**A. Dissolving the Solid**

1. Find a solvent with a steep solubility-vs.-temperature characteristic (done by trial and error using small amounts of material or by consulting a handbook).
2. Heat the desired solvent to its boiling point.
3. Dissolve the solid in a **minimum** of boiling solvent (either in a flask or a Craig tube).
4. If necessary, add decolorizing charcoal or decolorize the solution on a silica gel or alumina column.

B. Removing Insoluble Impurities

1. Decant or remove the solution with a Pasteur pipette, or
2. Filter the hot solution through a fluted filter, a filtering pipette, or a filter-tip pipette to remove insoluble impurities or charcoal.

NOTE: If no decolorizing charcoal has been added or if there are no undissolved particles, Part B should be omitted.

C. Crystallizing the Solid

1. Allow the solution to cool.
2. If crystals appear, cool the mixture in an ice-water bath (if desired) and go to Part D. If crystals do not appear, go to the next step.
3. Inducing crystallization
 - (a) Scratch the flask with a glass rod; or, if using a Craig tube, dip a glass rod or spatula into the solution, let the liquid evaporate, and place the glass rod or spatula back into the solution to seed it.
 - (b) Seed the solution with original solid, if available.
 - (c) Cool the solution in an ice-water bath.
 - (d) Evaporate excess solvent and allow the solution to cool again.

D. Collecting and Drying the Crystals

1. Collect crystals by vacuum filtration using a Hirsch funnel or by centrifugation using a Craig tube.
 2. If using a Hirsch funnel, rinse crystals with a small portion of **cold** solvent.
 3. Continue suction until crystals are nearly dry, if using vacuum filtration.
 4. Drying
 - (a) Air-dry the crystals, or
 - (b) Place the crystals in a drying oven, or
 - (c) Dry the crystals *in vacuo*.
-

PART D. ADDITIONAL EXPERIMENTAL CONSIDERATIONS: SEMIMICROSCALE AND MICROSCALE

11.5 Selecting a Solvent

A solvent that dissolves little of the material to be crystallized when it is cold but a great deal of the material when it is hot is a good solvent for crystallization. Quite often, correct crystallization solvents are indicated in the experimental procedures that you will be following. When a solvent is not specified in a procedure, you can determine a good crystallization solvent by consulting a handbook or making an educated guess based on polarities, both discussed in this section. A third approach, involving experimentation, is discussed in Section 11.6.

With compounds that are well known, the correct crystallization solvent has already been determined through the experiments of earlier researchers. In such

cases, the chemical literature can be consulted to determine which solvent should be used. Sources such as *The Merck Index* or the *CRC Handbook of Chemistry and Physics* may provide this information.


For example, consider naphthalene, which is found in *The Merck Index*, which states under the entry for naphthalene: “Monoclinic prismatic plates from ether.” This statement means that naphthalene can be crystallized from ether. It also gives the type of crystal structure. Unfortunately, the crystal structure may be given without reference to the solvent. Another way to determine the best solvent is by looking at solubility-vs.-temperature data. When this is given, a good solvent is one in which the solubility of the compound increases significantly as the temperature increases. Sometimes, the solubility data will be given for only cold solvent and boiling solvent. This should provide enough information to determine whether this would be a good solvent for crystallization.

In most cases, however, the handbooks will state only whether a compound is soluble or not in a given solvent, usually at room temperature. Determining a good solvent for crystallization from this information can be somewhat difficult. The solvent in which the compound is soluble may or may not be an appropriate solvent for crystallization. Sometimes, the compound may be too soluble in the solvent at all temperatures, and you would recover little of your product if this solvent were used for crystallization. It is possible that an appropriate solvent would be the one in which the compound is nearly insoluble at room temperature because the solubility-vs.-temperature curve is steep. Although the solubility information may give you some ideas about what solvents to try, you will most likely need to determine a good crystallizing solvent by experimentation as described in Section 11.6.

When using *The Merck Index* or *Handbook of Chemistry and Physics*, you should be aware that alcohol is frequently listed as a solvent. This generally refers to 95% or 100% ethyl alcohol. Because 100% (absolute) ethyl alcohol is more expensive than 95% ethyl alcohol, the cheaper grade is usually used in the chemistry laboratory. Another solvent frequently listed is benzene. Benzene is a known carcinogen, so it is rarely used in student laboratories. Toluene is a suitable substitute; the solubility behavior of a substance in benzene and toluene is so similar that you may assume any statement made about benzene also applies to toluene.

Another way to identify a solvent for crystallization is to consider the polarities of the compound and the solvents. Generally, you would look for a solvent that has a polarity somewhat similar to that of the compound to be crystallized. Consider the compound sulfanilamide, shown in the figure. There are several polar bonds in sulfanilamide, the NH and the SO bonds. In addition, the NH₂ groups and the oxygen atoms in sulfanilamide can form hydrogen bonds. Although the benzene ring portion of sulfanilamide is nonpolar, sulfanilamide has an intermediate polarity because of the polar groups. A common organic solvent of intermediate polarity is 95% ethyl alcohol. Therefore, it is likely that sulfanilamide would be soluble in 95% ethyl alcohol because they have similar polarities. (Note that the

Sulfanilamide

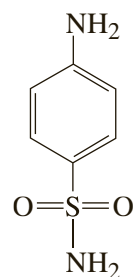

other 5% in 95% ethyl alcohol is usually a substance such as water or isopropyl alcohol, which does not alter the overall polarity of the solvent.) Although this kind of analysis is a good first step in determining an appropriate solvent for crystallization, without more information it is not enough to predict the shape of the solubility curve for the temperature-vs.-solubility data (see Figure 11.1). Therefore, knowing that sulfanilamide is soluble in 95% ethyl alcohol does not necessarily mean that this is a good solvent for crystallizing sulfanilamide. You would still need to test the solvent to see if it is appropriate. The solubility curve for sulfanilamide (see Figure 11.2) indicates that 95% ethyl alcohol is a good solvent for crystallizing this substance.

When choosing a crystallization solvent, do not select one whose boiling point is higher than the melting point of the substance (solute) to be crystallized. If the boiling point of the solvent is too high, the substance may come out of solution as a liquid rather than a crystalline solid. In such a case, the solid may **oil out**. Oiling out occurs when, on cooling the solution to induce crystallization, the solute begins to come out of solution at a temperature above its melting point. The solute will then come out of solution as a liquid. Furthermore, as cooling continues, the substance may still not crystallize; rather, it will become a supercooled liquid. Oils may eventually solidify if the temperature is lowered, but often they will not actually crystallize. Instead, the solidified oil will be an amorphous solid or a hardened mass. In this case, purification of the substance will not have occurred as it does when the solid is crystalline. It can be difficult to deal with oils when trying to obtain a pure substance. You must try to redissolve them and hope that the substance will crystallize with slow, careful cooling. During the cooling period, it may be helpful to scratch the glass container where the oil is present with a glass stirring rod that has not been fire-polished. Seeding the oil as it cools with a small sample of the original solid is another technique that is sometimes helpful in working with difficult oils. Other methods of inducing crystallization are discussed in Section 11.8.

One additional criterion for selecting the correct crystallization solvent is the **volatility** of that solvent. Volatile solvents have low boiling points or evaporate easily. A solvent with a low boiling point may be removed from the crystals through evaporation without much difficulty. It will be difficult to remove a solvent with a high boiling point from the crystals without heating them under vacuum. On the other hand, solvents with very low boiling points are not ideal for crystallizations. The recovery will not be as great with low boiling solvents because they cannot be heated past the boiling point and the amount of cooling that occurs during the crystallization step will be relatively small. Diethyl ether (bp = 35°C) and methylene chloride (bp = 41°C) are not often used as crystallization solvents.

Table 11.2 lists common crystallization solvents. The solvents used most commonly are listed in the table first.

11.6 Testing Solvents for Crystallization

When the appropriate solvent is not known, select a solvent for crystallization by experimenting with various solvents and a small amount of the material to be crystallized. Experiments are conducted on a small test tube scale before the entire quantity of material is committed to a particular solvent. Such trial-and-error methods are common when trying to purify a solid material that has not been studied.

Procedure

1. Place about 0.05 g of the sample in a test tube.
2. Add about 0.5 mL of solvent at room temperature and stir the mixture by rapidly twirling a microspatula between your fingers. If all (or almost all) of the

Table 11.2 Common solvents for crystallization

	Boils (°C)	Freezes (°C)	Soluble in H ₂ O	Flammability
Water	100	0	+	–
Methanol	65	*	+	+
95% Ethanol	78	*	+	+
Ligroin	60–90	*	–	+
Toluene	111	*	–	+
Chloroform**	61	*	–	–
Acetic acid	118	17	+	+
Dioxane**	101	11	+	+
Acetone	56	*	+	+
Diethyl ether	35	*	Slightly	++
Petroleum ether	30–60	*	–	++
Methylene chloride	41	*	–	–
Carbon tetrachloride**	77	*	–	–

*Lower than 0°C (ice temperature).

**Suspected carcinogen.

solid dissolves at room temperature, then your solid is *probably* too soluble in this solvent and little compound would be recovered if this solvent were used. Select another solvent.

3. If none (or very little) of the solid dissolves at room temperature, heat the tube carefully and stir with a spatula. (A hot water bath is perhaps better than an aluminum block because you can more easily control the temperature of the hot water bath. The temperature of the hot water bath should be slightly higher than the boiling point of the solvent.) Add more solvent dropwise while continuing to heat and stir. Continue adding solvent until the solid dissolves, but do not add more than about 1.5 mL (total) of solvent. If all the solid dissolves, go to step 4. If all the solid has not dissolved by the time you have added 1.5 mL of solvent, this is probably not a good solvent. However, if most of the solid has dissolved at this point, you might try adding a little more solvent. Remember to heat and stir at all times during this step.
4. If the solid dissolves in about 1.5 mL or less of boiling solvent, then remove the test tube from the heat source, stopper the tube, and allow it to cool to room temperature. Then place it in an ice-water bath. If a lot of crystals come out, this is most likely a good solvent. If crystals do not come out, scratch the sides of the tube with a glass stirring rod to induce crystallization. If crystals still do not form, this is probably not a good solvent.

Comments about This Procedure

1. Selecting a good solvent is something of an art. There is no perfect procedure that can be used in all cases. You must think about what you are doing and use some common sense in deciding whether to use a particular solvent.
2. Do not heat the mixture above the melting point of your solid. This can occur most easily when the boiling point of the solvent is higher than the melting

point of the solid. Normally, do not select a solvent that has a higher boiling point than the melting point of the substance. If you do, make certain that you do not heat the mixture beyond the melting point of your solid.

11.7 Decolorization

Small amounts of highly colored impurities may make the original crystallization solution appear colored; this color can often be removed by **decolorization**, either by using activated charcoal (often called Norit) or by passing the solution through a column packed with alumina or silica gel. A decolorizing step should be performed only if the color is due to impurities, not to the color of the desired product, and if the color is significant. Small amounts of colored impurities will remain in solution during crystallization, making the decolorizing step unnecessary. The use of activated charcoal is described separately for macroscale and microscale crystallizations, and then the column technique, which can be used with both crystallization techniques, is described.

A. Semimicroscale—Powdered Charcoal

As soon as the solute is dissolved in the minimum amount of boiling solvent, the solution is allowed to cool slightly, and a small amount of Norit (powdered charcoal) is added to the mixture. The Norit adsorbs the impurities. When performing a crystallization in which the filtration is performed with a fluted filter, you should add powdered Norit because it has a larger surface area and can remove impurities more effectively. A reasonable amount of Norit is what could be held on the end of a microspatula, or about 0.01–0.02 g. If too much Norit is used, it will adsorb product as well as impurities. A small amount of Norit should be used, and its use should be repeated if necessary. (It is difficult to determine if the initial amount added is sufficient until after the solution is filtered because the suspended particles of charcoal will obscure the color of the liquid.) Caution should be exercised so that the solution does not froth or erupt when the finely divided charcoal is added. The mixture is boiled with the Norit for several minutes and then filtered by gravity, using a fluted filter (see Section 11.3 and Technique 8, Section 8.1B), and the crystallization is carried forward as described in Section 11.3.

The Norit preferentially adsorbs the colored impurities and removes them from the solution. The technique seems to be most effective with hydroxylic solvents. In using Norit, be careful not to breathe the dust. Normally, small quantities are used so that little risk of lung irritation exists.

B. Microscale—Pelletized Norit

If the crystallization is being performed in a Craig tube, it is advisable to use pelletized Norit. Although this is not as effective in removing impurities as powdered Norit, it is easier to remove, and the amount of pelletized Norit required is more easily determined because you can see the solution as it is being decolorized. Again, the Norit is added to the hot solution (the solution should not be boiling) after the solid has dissolved. This should be performed in a test tube rather than in a Craig tube. About 0.02 g is added, and the mixture is boiled for a minute or so to see if more Norit is required. More Norit is added, if necessary, and the liquid is boiled again. It is important not to add too much pelletized Norit because the Norit will also adsorb some of the desired material, and it is possible that not all the color can be removed no matter how much is added. The decolorized solution is then removed with a preheated filter-tip pipette (see Technique 8, Section 8.6) to filter the mixture and transferred to a Craig tube for crystallization as described in Section 11.4.

C. Decolorization on a Column

The other method for decolorizing a solution is to pass the solution through a column containing alumina or silica gel. The adsorbent removes the colored impurities while allowing the desired material to pass through (see Technique 8, Figure 8.6, and Technique 19, Section 19.14). If this technique is used, it will be necessary to dilute the solution with additional solvent to prevent crystallization from occurring during the process. The excess solvent must be evaporated after the solution is passed through the column (Technique 7, Section 7.10), and the crystallization procedure is continued as described in Sections 11.3 or 11.4.

11.8 Inducing Crystallization

If a cooled solution does not crystallize, several techniques may be used to induce crystallization. Although identical in principle, the actual procedures vary slightly when macroscale and microscale crystallizations are performed.

A. Semimicroscale

In the first technique, you should try scratching the inside surface of the flask vigorously with a glass rod that has not been fire-polished. The motion of the rod should be vertical (in and out of the solution) and should be vigorous enough to produce an audible scratching. Such scratching often induces crystallization, although the effect is not well understood. The high-frequency vibrations may have something to do with initiating crystallization; or perhaps—a more likely possibility—small amounts of solution dry by evaporation on the side of the flask, and the dried solute is pushed into the solution. These small amounts of material would provide “seed crystals,” or nuclei, on which crystallization may begin.

A second technique that can be used to induce crystallization is to cool the solution in an ice bath. This method decreases the solubility of the solute.

A third technique is useful when small amounts of the original material to be crystallized are saved. The saved material can be used to “seed” the cooled solution. A small crystal dropped into the cooled flask often will start the crystallization—this is called **seeding**.

If all these measures fail to induce crystallization, it is likely that too much solvent was added. The excess solvent must then be evaporated (Technique 7, Section 7.10) and the solution allowed to cool.

B. Microscale

The strategy is basically the same as described for macroscale crystallizations. Scratching vigorously with a glass rod should be avoided, however, because the Craig tube is fragile and expensive. Scratching *gently* is allowed.

Another measure is to dip a spatula or glass stirring rod into the solution and allow the solvent to evaporate so that a small amount of solid will form on the surface of the spatula or glass rod. When placed back into the solution, the solid will seed the solution. A small amount of the original material, if some was saved, may also be used to seed the solution.

A third technique is to cool the Craig tube in an ice-water bath. This method may also be combined with either of the previous suggestions.

If none of these measures is successful, it is possible that too much solvent is present, and it may be necessary to evaporate some of the solvent (Technique 7, Section 7.10) and allow the solution to cool again.

11.9 Drying Crystals

The most common method of drying crystals involves allowing them to dry in air. Several methods are illustrated in Figure 11.7, below. In all three methods, the crystals must be covered to prevent accumulation of dust particles. Note that in each method, the spout on the beaker provides an opening so that solvent vapor can escape from the system. The advantage of this method is that heat is not required, thus reducing the danger of decomposition or melting; however, exposure to atmospheric moisture may cause the hydration of strongly hygroscopic materials. A **hygroscopic** substance is a substance that absorbs moisture from the air.

Another method of drying crystals is to place the crystals on a watch glass, a clay plate, or a piece of absorbent paper in an oven. Although this method is simple, some possible difficulties deserve mention. Crystals that sublime readily should not be dried in an oven because they might vaporize and disappear. Care should be taken that the temperature of the oven does not exceed the melting point of the crystals. Remember that the melting point of crystals is lowered by the presence of solvent; allow for this melting-point depression when selecting a suitable oven temperature. Some materials decompose on exposure to heat, and they should not be dried in an oven. Finally, when many different samples are being dried in the same oven, crystals might be lost due to confusion or reaction with another person's sample. It is important to label the crystals when they are placed in the oven.

A third method, which requires neither heat nor exposure to atmospheric moisture, is drying *in vacuo*. Two procedures are illustrated in Figure 11.8.

Procedure A

In this method, a desiccator is used. The sample is placed under vacuum in the presence of a drying agent. Two potential problems must be noted. The first deals with samples that sublime readily. Under vacuum, the likelihood of sublimation is increased. The second problem deals with the vacuum desiccator itself. Because the surface area of glass that is under vacuum is large, there is some danger that the desiccator could implode. A vacuum desiccator should never be used unless it has been placed within a protective metal container (cage). If a cage is not available, the desiccator can be wrapped with electrical or duct tape. If you use an aspirator as a source of vacuum, you should use a water trap (see Technique 8, Figure 8.5).

Procedure B

This method can be accomplished with a round-bottom flask and a thermometer adapter equipped with a short piece of glass tubing, as illustrated in Figure 11.8B.

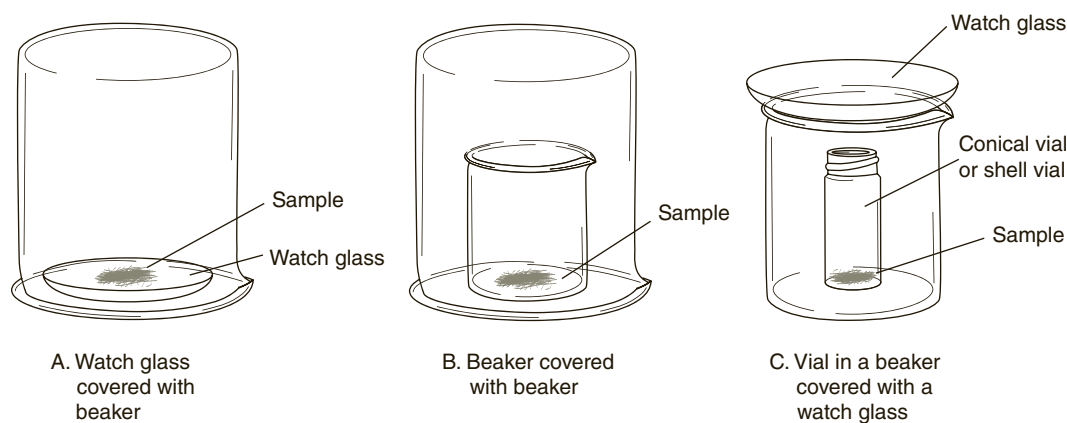


Figure 11.7
Methods for drying crystals in air.

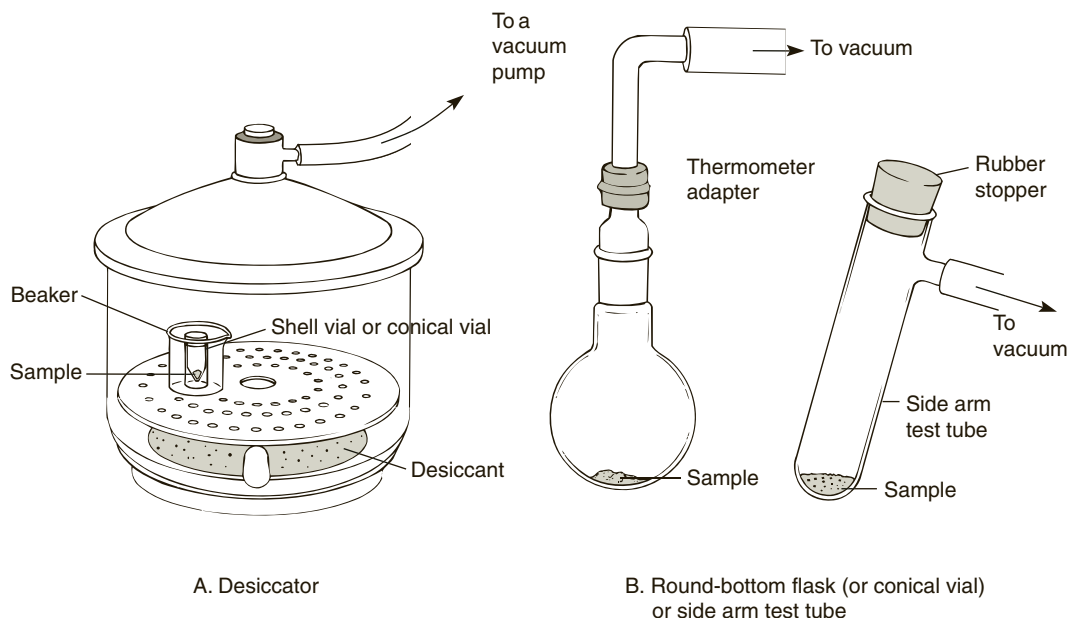


Figure 11.8
Methods for drying crystals in a vacuum.

In microscale work, the apparatus with the round-bottom flask can be modified by replacing the round-bottom flask with a conical vial. The glass tubing is connected by vacuum tubing to either an aspirator or a vacuum pump. A convenient alternative, using a side arm test tube, is also shown in Figure 11.8B. With either apparatus, install a water trap when an aspirator is used.

11.10 Mixed Solvents

Often, the desired solubility characteristics for a particular compound are not found in a single solvent. In these cases, a mixed solvent may be used. You simply select a first solvent in which the solute is soluble and a second solvent, miscible with the first, in which the solute is relatively insoluble. The compound is dissolved in a minimum amount of the boiling solvent in which it is soluble. Following this, the second hot solvent is added to the boiling mixture, dropwise, until the mixture barely becomes cloudy. The cloudiness indicates precipitation. At this point, more of the first solvent should be added. Just enough is added to clear the cloudy mixture. At that point, the solution is saturated, and as it cools, crystals should separate. Common solvent mixtures are listed in Table 11.3.

It is important not to add an excess of the second solvent or to cool the solution too rapidly. Either of these actions may cause the solute to oil out, or separate as a viscous liquid. If this happens, reheat the solution and add more of the first solvent.

Table 11.3 Common solvent pairs for crystallization

Methanol–water	Ether–acetone
Ethanol–water	Ether–petroleum ether
Acetic acid–water	Toluene–ligroin
Acetone–water	Methylene chloride–methanol
Ether–methanol	Dioxane ^a –water

^aSuspected carcinogen.

PROBLEMS

1. Listed below are solubility-vs.-temperature data for an organic substance A dissolved in water.

Temperature (°C)	Solubility of A in 100 mL of Water (g)
0	1.5
20	3.0
40	6.5
60	11.0
80	17.0

- Graph the solubility of A vs. temperature. Use the data given in the table. Connect the data points with a smooth curve.
 - Suppose 0.1 g of A and 1.0 mL of water were mixed and heated to 80°C. Would all the substance A dissolve?
 - The solution prepared in (b) is cooled. At what temperature will crystals of A appear?
 - Suppose the cooling described in (c) were continued to 0°C. How many grams of A would come out of solution? Explain how you obtained your answer.
- What would likely happen if a hot saturated solution were filtered by vacuum filtration using a Hirsch funnel? (*Hint*: The mixture will cool as it comes in contact with the Hirsch funnel.)
 - A compound you have prepared is reported in the literature to have a pale yellow color. When the substance is dissolved in hot solvent to purify it by crystallization, the resulting solution is yellow. Should you use decolorizing charcoal before allowing the hot solution to cool? Explain your answer.
 - After a crude product is dissolved in 1.5 mL of hot solvent, the resulting solution is dark brown. Because the pure compound is reported in the literature to be colorless, it is necessary to perform a decolorizing procedure. Should you use pelletized Norit or powdered activated charcoal to decolorize the solution? Explain your answer.
 - While performing a crystallization, you obtain a light tan solution after dissolving your crude product in hot solvent. A decolorizing step is determined to be unnecessary, and there are no solid impurities present. Should you perform a filtration to remove impurities before allowing the solution to cool? Why or why not?
 - Draw a graph of a cooling curve (temperature vs. time) for a solution of a solid substance that shows no supercooling effects. Assume that the solvent does not freeze.
 - Repeat the instructions in (a) for a solution of a solid substance that shows some supercooling behavior but eventually yields crystals if the solution is cooled sufficiently.
 - A solid substance A is soluble in water to the extent of 10 mg/mL of water at 25°C and 100 mg/mL of water at 100°C. You have a sample that contains 100 mg of A and an impurity B. See Section 11.2 for guidance in answering these questions.

- a. Assuming that 2 mg of B are present along with 100 mg of A, describe how you can purify A if B is completely insoluble in water. Your description should include the volume of solvent required.
 - b. Assuming that 2 mg of the impurity B are present along with 100 mg of A, describe how you can purify A if B has the same solubility behavior as A. Will one crystallization produce pure A? (Assume that the solubilities of both A and B are unaffected by the presence of the other substance.)
 - c. Assume that 25 mg of the impurity B are present along with 100 mg of A. Describe how you can purify A if B has the same solubility behavior as A. Each time, use the minimum amount of water to just dissolve the solid. Will one crystallization produce absolutely pure A? How many crystallizations would be needed to produce pure A? How much A will have been recovered when the crystallizations have been completed?
8. Consider the crystallization of sulfanilamide from 95% ethyl alcohol. If impure sulfanilamide is dissolved in the minimum amount of 95% ethyl alcohol at 40°C rather than 78°C (the boiling point of ethyl alcohol), how would this affect the percentage of recovery of pure sulfanilamide? Explain your answer.

Extractions, Separations, and Drying Agents

PART A. THEORY

12.1 Extraction

Transferring a solute from one solvent into another is called **extraction**, or, more precisely, liquid–liquid extraction. The solute is extracted from one solvent into the other because the solute is more soluble in the second solvent than in the first. The two solvents must not be **miscible** (mix freely), and they must form two separate **phases**, or layers, in order for this procedure to work. Extraction is used in many ways in organic chemistry. Many natural products (organic chemicals that exist in nature) are present in animal and plant tissues having high water content. Extracting these tissues with a water-immiscible solvent is useful for isolating natural products. Often, diethyl ether (commonly referred to as “ether”) is used for this purpose. Sometimes alternative water-immiscible solvents such as hexane, petroleum ether, ligroin, and methylene chloride are used. For instance, caffeine, a natural product, can be extracted from an aqueous tea solution by shaking it successively with several portions of methylene chloride.

A generalized extraction process that uses a conical vial is illustrated in Figure 12.1. The first solvent contains a mixture of black and white molecules (Figure 12.1A). A second solvent that is not miscible with the first is added.

After the vial is capped and shaken, the layers separate. In this example, the second solvent is less dense, so it becomes the top layer (Figure 12.1B). Because of differences in physical properties, the white molecules are more soluble in the second solvent, whereas the black molecules are more soluble in the first solvent. Most of the white molecules are in the upper layer, but there are some black molecules there, too. Likewise, most of the black molecules are in the lower layer. However, there are a few white molecules in this lower phase. A Pasteur pipette may be used to remove the lower layer (Figure 12.1C). In this way, a partial separation of black and white molecules has been achieved. In this example, notice that it was not possible to effect a complete separation with one extraction. This is a common occurrence in organic chemistry. Many organic substances are soluble in both water and organic solvents.

Water can be used to extract or “wash” water-soluble impurities from an organic reaction mixture. To carry out a “washing” operation, you add water to the reaction mixture contained in a conical vial. After capping the vial and shaking it, you allow the organic layer and the aqueous (water) layer to separate from each other in the vial. A water wash removes highly polar and water-soluble materials, such as sulfuric acid, hydrochloric acid, or sodium hydroxide from the organic layer. A water wash can also be used to remove water-soluble and low-molecular-weight compounds, such as ethanol or acetic acid, from the

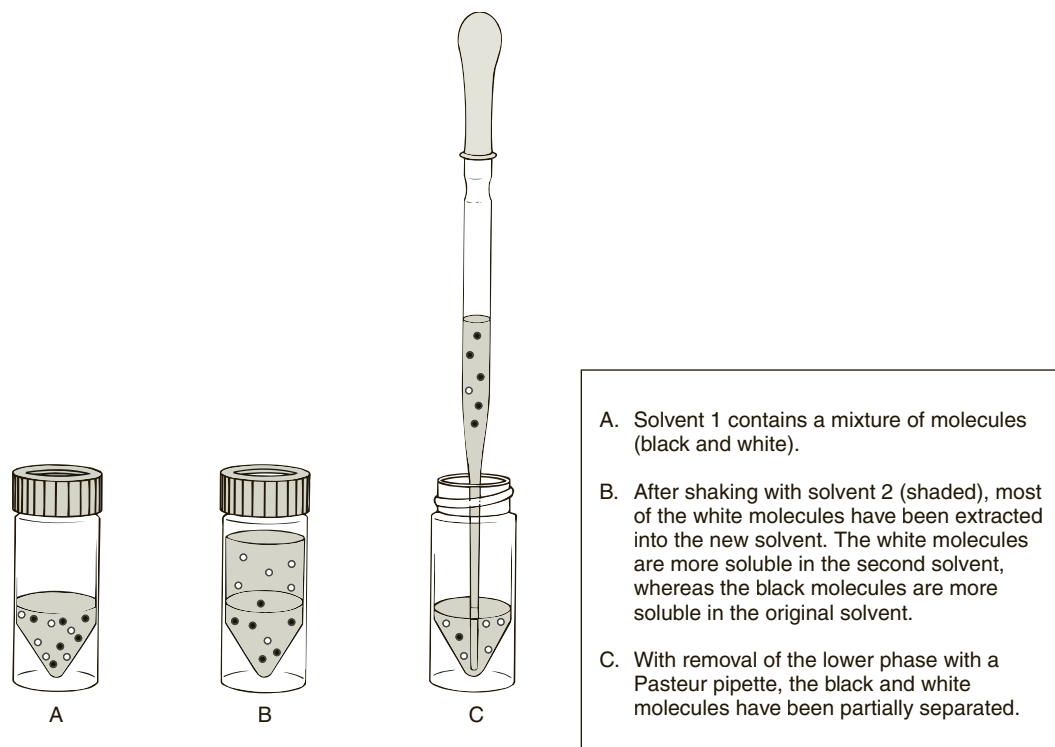


Figure 12.1
The extraction process.

organic layer. The washing operation helps purify the desired organic compound present in the original reaction mixture.

12.2 Distribution Coefficient

When a solution (solute A in solvent 1) is shaken with a second solvent (solvent 2) with which it is not miscible, the solute distributes itself between the two liquid phases. When the two phases have separated again into two distinct solvent layers, an equilibrium will have been achieved such that the ratio of the concentrations of the solute in each layer defines a constant. The constant, called the **distribution coefficient** (or partition coefficient) K , is defined by

$$K = \frac{C_2}{C_1}$$

where C_1 and C_2 are the concentrations at equilibrium, in grams per liter or milligrams per milliliter of solute A in solvent 1 and in solvent 2, respectively. This relationship is a ratio of two concentrations and is independent of the actual amounts of the two solvents mixed. The distribution coefficient has a constant value for each solute considered and depends on the nature of the solvents used in each case.

Not all the solute will be transferred to solvent 2 in a single extraction unless K is very large. Usually it takes several extractions to remove all the solute from solvent 1. In extracting a solute from a solution, it is always better to use several small portions of the second solvent than to make a single extraction with a large portion. Suppose that, as an illustration, a particular extraction proceeds with a distribution coefficient of 10. The system consists of 50 mg of organic compound dissolved in 1.00 mL of water (solvent 1). In this illustration, the effectiveness of three 0.50-mL extractions with ether (solvent 2) is compared with one 1.50-mL extraction with ether. In the first 0.50-mL extraction, the amount extracted into the ether layer is

given by the following calculation. The amount of compound remaining in the aqueous phase is given by x .

$$K = 10 = \frac{C_2}{C_1} = \frac{\left(\frac{(50.0 - x) \text{ mg}}{0.50 \text{ mL ether}}\right)}{\left(\frac{x \text{ mg}}{1.00 \text{ mL water}}\right)}; 10 = \frac{(50.0 - x)(1.00)}{0.50x}$$

$$5.0x = 50.0 - x$$

$$6.0x = 50.0$$

$$x = 8.3 \text{ mg remaining in the aqueous layer}$$

$$50.0 - x = 41.7 \text{ mg in the ether layer}$$

As a check on the calculation, it is possible to substitute the value 8.3 mg for x in the original equation and demonstrate that the concentration in the ether phase divided by the concentration in the water phase equals the distribution coefficient.

$$\frac{\left(\frac{(50.0 - x) \text{ mg}}{0.50 \text{ mL ether}}\right) \frac{41.7}{0.50}}{\left(\frac{x \text{ mg}}{1.00 \text{ mL water}}\right) \frac{8.3}{1.00}} = \frac{83 \text{ mg/mL}}{8.3 \text{ mg/mL}} = 10 = K$$

The second extraction with another 0.50-mL portion of fresh ether is performed on the aqueous phase, which now contains 8.3 mg of the solute. The amount of solute extracted is given by the calculation shown in Figure 12.2. Also shown in the



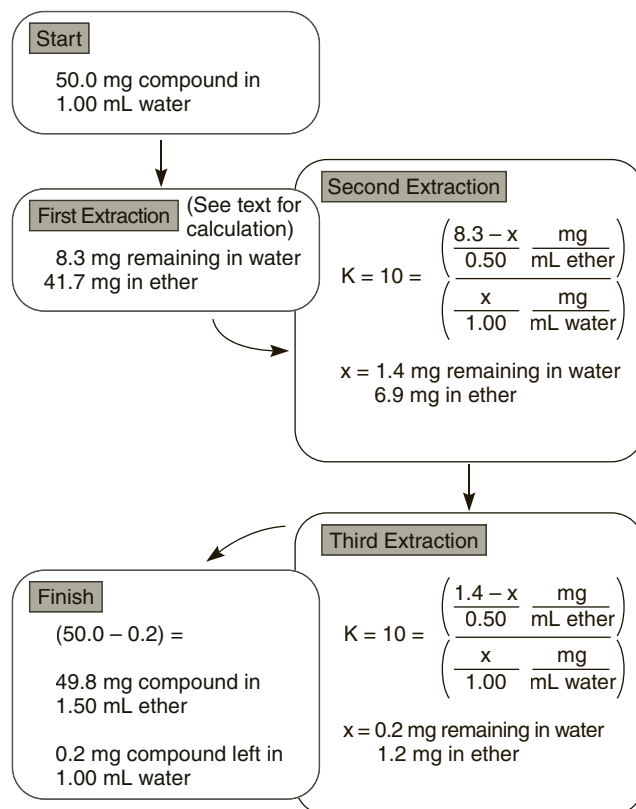


Figure 12.2

The result of extraction of 50.0 mg of compound in 1.00 mL of water by three successive 0.50-mL portions of ether. Compare this result with that of Figure 12.3.

**Figure 12.3**

The result of extraction of 50.0 mg of compound in 1.00 mL of water with one 1.5-mL portion of ether. Compare this result with that of Figure 12.2.

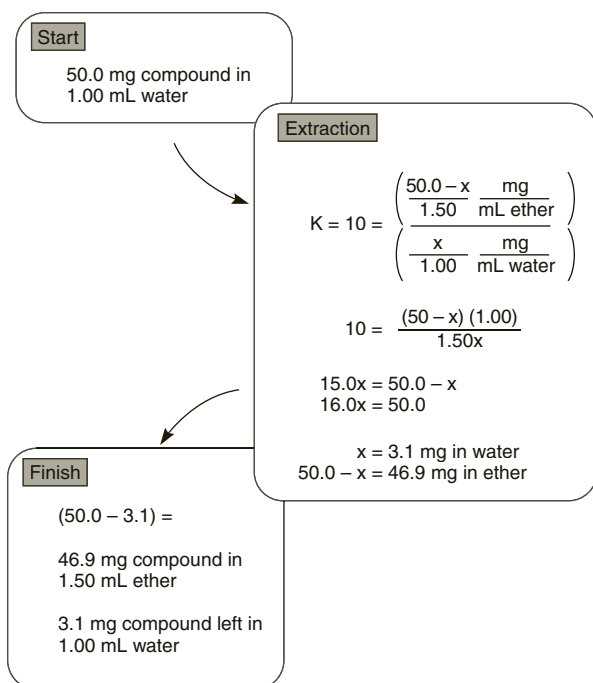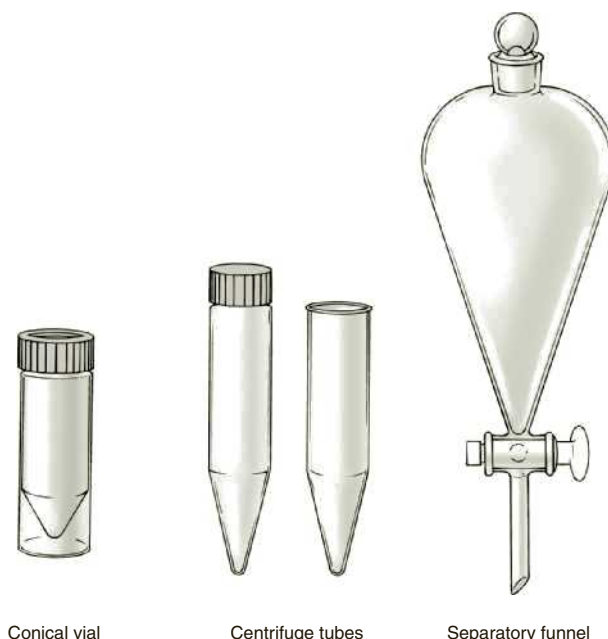

Conical vial

Centrifuge tubes

Separatory funnel

Figure 12.4

Apparatus used for extraction.

figure is a calculation for a third extraction with another 0.50-mL portion of ether. This third extraction will transfer 1.2 mg of solute into the ether layer, leaving 0.2 mg of solute remaining in the water layer. A total of 49.8 mg of solute will be extracted into the combined ether layers, and 0.2 mg will remain in the aqueous phase.

Figure 12.3 shows the result of a *single* extraction with 1.50 mL of ether. As shown there, 46.9 mg of solute was extracted into the ether layer, leaving 3.1 mg of compound in the aqueous phase. One can see that three successive 0.50-mL ether extractions (Figure 12.2) succeeded in removing 2.9 mg more solute from the aqueous phase than using one 1.50-mL portion of ether (Figure 12.3). This differential represents 5.8% of the total material.

12.3 Choosing an Extraction Method and a Solvent

Three types of apparatus are used for extractions: conical vials, centrifuge tubes, and separatory funnels. These are shown in Figure 12.4. Conical vials may be used with volumes of less than 4 mL; volumes of up to 10 mL may be handled in centrifuge tubes. A centrifuge tube equipped with a screw cap is particularly useful for extractions. The separatory funnel is used in large-scale reactions. Each type of equipment is discussed in a separate section.

Before using a conical vial for an extraction, make sure that the capped conical vial does not leak when shaken. To do this, place some water in the conical vial, place the Teflon liner in the cap, and screw the cap securely onto the conical vial. Shake the vial vigorously and check for leaks. Conical vials that are used for extractions must not be chipped on the edge of the vial or they will not seal adequately. If there is a leak, try tightening the cap or replacing the Teflon liner with another one. Sometimes it helps to use the silicone rubber side of the liner to seal the conical vial. Some laboratories are supplied with Teflon stoppers that fit into the 5-mL conical vials. You may find that this stopper eliminates leakage.

Table 12.1 Densities of common extraction solvents

<i>Solvent</i>	<i>Density (g/mL)</i>
Ligroin	0.67–0.69
Diethyl ether	0.71
Toluene	0.87
Water	1.00
Methylene chloride	1.330

When shaking the conical vial, do it gently at first in a rocking motion. When it is clear that an emulsion will not form (see Section 12.10), you can shake it more vigorously.

In some cases, adequate mixing can be achieved by spinning your microspatula for at least 10 minutes in the conical vial. Another technique of mixing involves drawing up the mixture into a Pasteur pipette and squirting it rapidly back into the vial. Repeat this process for at least 5 minutes to obtain an adequate extraction.

If you are using a screw-cap centrifuge tube, put some water in the tube, cap it, and shake it vigorously to check for leaks. If the centrifuge tube leaks, try replacing the cap with another one. If available in the laboratory, a vortex mixer may be used to mix the phases. A vortex mixer works well with a variety of containers, including small flasks, test tubes, conical vials, and centrifuge tubes. You start the mixing action on a vortex mixer by holding the test tube or other container on one of the pads. The unit mixes the sample by high-frequency vibration.

Most extractions consist of an aqueous phase and an organic phase. In order to extract a substance from an aqueous phase, an organic solvent that is not miscible with water must be used. Table 12.1 lists a number of the common organic solvents that are not miscible with water and are used for extraction.

Those solvents that have a density less than that of water (1.00 g/mL) will separate as the top layer when shaken with water. Those solvents that have a greater density than water will separate into the lower layer. For instance, diethyl ether ($d = 0.71$ g/mL) when shaken with water will form the upper layer, whereas methylene chloride ($d = 1.33$ g/mL) will form the lower layer. When an extraction is performed, slightly different methods are used when you wish to separate the lower layer (whether it is the aqueous layer or the organic layer) than when you wish to separate the upper layer.


PART B. MICROSCALE EXTRACTION

12.4 The Conical Vial—Separating the Lower Layer

The 5-mL conical vial is the most useful piece of equipment for carrying out extractions on a microscale level. In this section, we consider the method for removing the lower layer. A concrete example would be the extraction of a desired product from an aqueous layer using methylene chloride ($d = 1.33$ g/mL) as the extraction solvent. Methods for removal of the upper layer are discussed in the next section.

NOTE: Always place a conical vial in a small beaker to prevent the vial from falling over.

Removing the Lower Layer. Suppose that we extract an aqueous solution with methylene chloride. This solvent is denser than water and will settle to the bottom

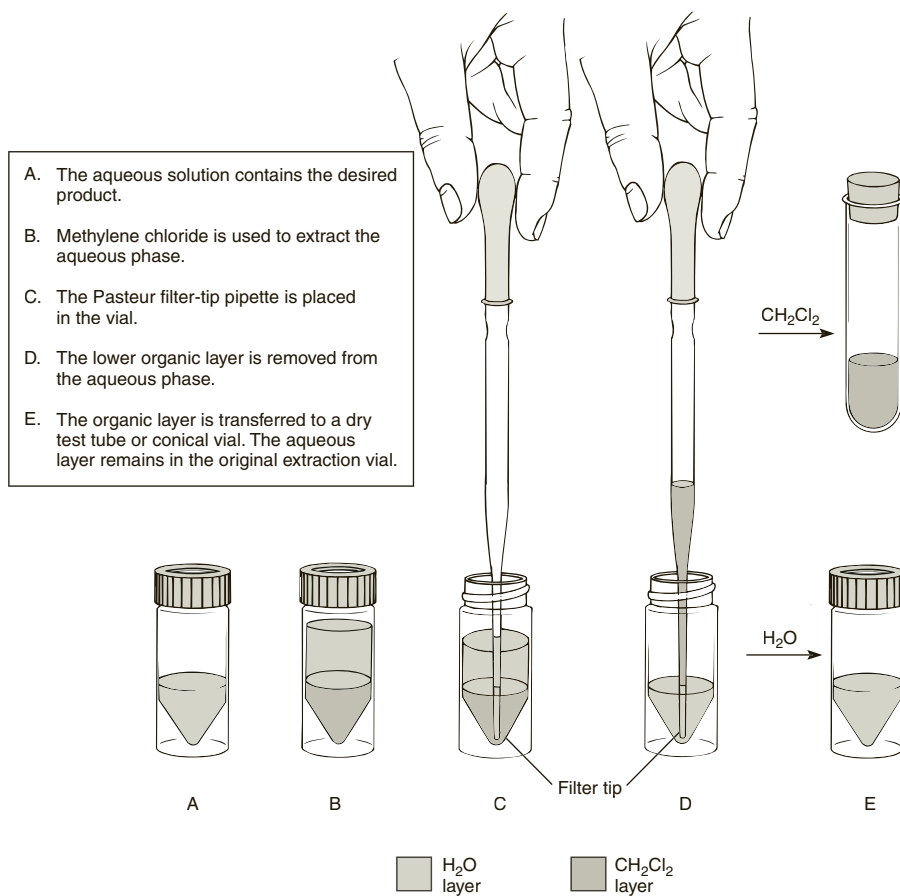**Figure 12.5**

Extraction of an aqueous solution using a solvent denser than water: methylene chloride.

of the conical vial. Use the following procedure, which is illustrated in Figure 12.5, to remove the lower layer.

1. Place the aqueous phase containing the dissolved product into a 5-mL conical vial (Figure 12.5A).
2. Add about 1 mL of methylene chloride, cap the vial, and shake the mixture gently at first in a rocking motion and then more vigorously when it is clear that an emulsion will not form. Vent or unscrew the cap slightly to release the pressure in the vial. Allow the phases to separate completely so that you can detect two distinct layers in the vial. The organic phase will be the lower layer in the vial (Figure 12.5B). If necessary, tap the vial with your finger or stir the mixture gently if some of the organic phase is suspended in the aqueous layer.
3. Prepare a Pasteur filter-tip pipette (Technique 8, Section 8.6) using a 5 $\frac{3}{4}$ -inch pipette. Attach a 2-mL rubber bulb to the pipette, depress the bulb, and insert the pipette into the vial so that the tip touches the bottom (Figure 12.5C). The filter-tip pipette gives you better control in removing the lower layer. In some cases, however, you may be able to use a Pasteur pipette (no filter tip), but considerably more care must be taken to avoid losing liquid from the pipette during the transfer operation. With experience, you should be able to judge how much to squeeze the bulb to draw in the desired volume of liquid.
4. Slowly draw the lower layer (methylene chloride) into the pipette in such a way that you exclude the aqueous layer and any emulsion (Section 12.10) that

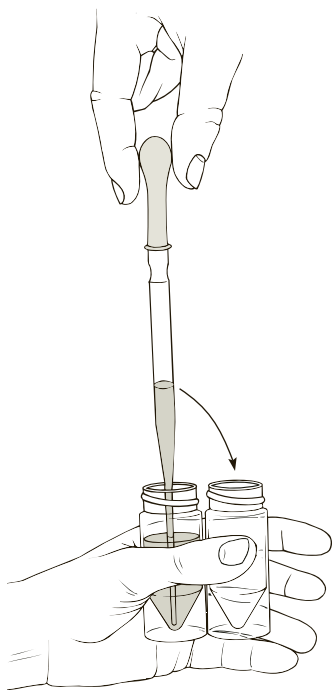


Figure 12.6
Method of holding vials
while transferring liquids.

might be at the interface between the layers (Figure 12.5D). Be sure to keep the tip of the pipette squarely in the V at the bottom of the vial.

5. Transfer the withdrawn organic phase into a *dry* test tube or another *dry* conical vial if one is available. It is best to have the test tube or vial located next to the extraction vial. Hold the vials in the same hand between your index finger and thumb, as shown in Figure 12.6. This avoids messy and disastrous transfers. The aqueous layer (upper layer) is left in the original conical vial (Figure 12.5E).

In performing an actual extraction in the laboratory, you would extract the aqueous phase with a second 1-mL portion of fresh methylene chloride to achieve a more complete extraction. Steps 2–5 would be repeated, and the organic layers from both extractions would be combined. In some cases, you may need to extract a third time with yet another 1-mL portion of methylene chloride. Again, the methylene chloride would be combined with the other extracts. The overall process would use three 1-mL portions of methylene chloride to transfer the product from the water layer into methylene chloride. Sometimes you will see the statement “extract the aqueous phase with three 1-mL portions of methylene chloride” in an experimental procedure. This statement describes in a shorter fashion the process described previously. Finally, the methylene chloride extracts will contain some water and must be dried with a drying agent as indicated in Section 12.9.

In this example, we extracted water with the heavy solvent methylene chloride and removed it as the lower layer. If you were extracting a light solvent (for instance, diethyl ether) with water, and you wished to keep the water layer, the water would be the lower layer and would be removed using the same procedure. You would not dry the water layer, however.


12.5 The Conical Vial—Separating the Upper Layer

In this section, we consider the method used when you wish to remove the upper layer. A concrete example would be the extraction of a desired product from an aqueous layer using diethyl ether ($d = 0.71$ g/mL) as the extraction solvent. Methods for removing the lower layer were discussed previously.

NOTE: Always place a conical vial in a small beaker to prevent the vial from falling over.

Removing the Upper Layer. Suppose we extract an aqueous solution with diethyl ether (ether). This solvent is less dense than water and will rise to the top of the conical vial. Use the following procedure, which is illustrated in Figure 12.7, to remove the upper layer.

1. Place the aqueous phase containing the dissolved product in a 5-mL conical vial (Figure 12.7A).
2. Add about 1 mL of ether, cap the vial, and shake the mixture vigorously. Vent or unscrew the cap slightly to release the pressure in the vial. Allow the phases to separate completely so that you can detect two distinct layers in the vial. The ether phase will be the upper layer in the vial (Figure 12.7B).
3. Prepare a Pasteur filter-tip pipette (Technique 8, Section 8.6) using a 5 $\frac{3}{4}$ -inch pipette. Attach a 2-mL rubber bulb to the pipette, depress the bulb, and insert the pipette into the vial so that the tip touches the bottom. The filter-tip pipette gives you better control in removing the lower layer. In some cases, however, you may be able to use a Pasteur pipette (no filter tip), but considerably more care must be taken to avoid losing liquid from the pipette during the transfer

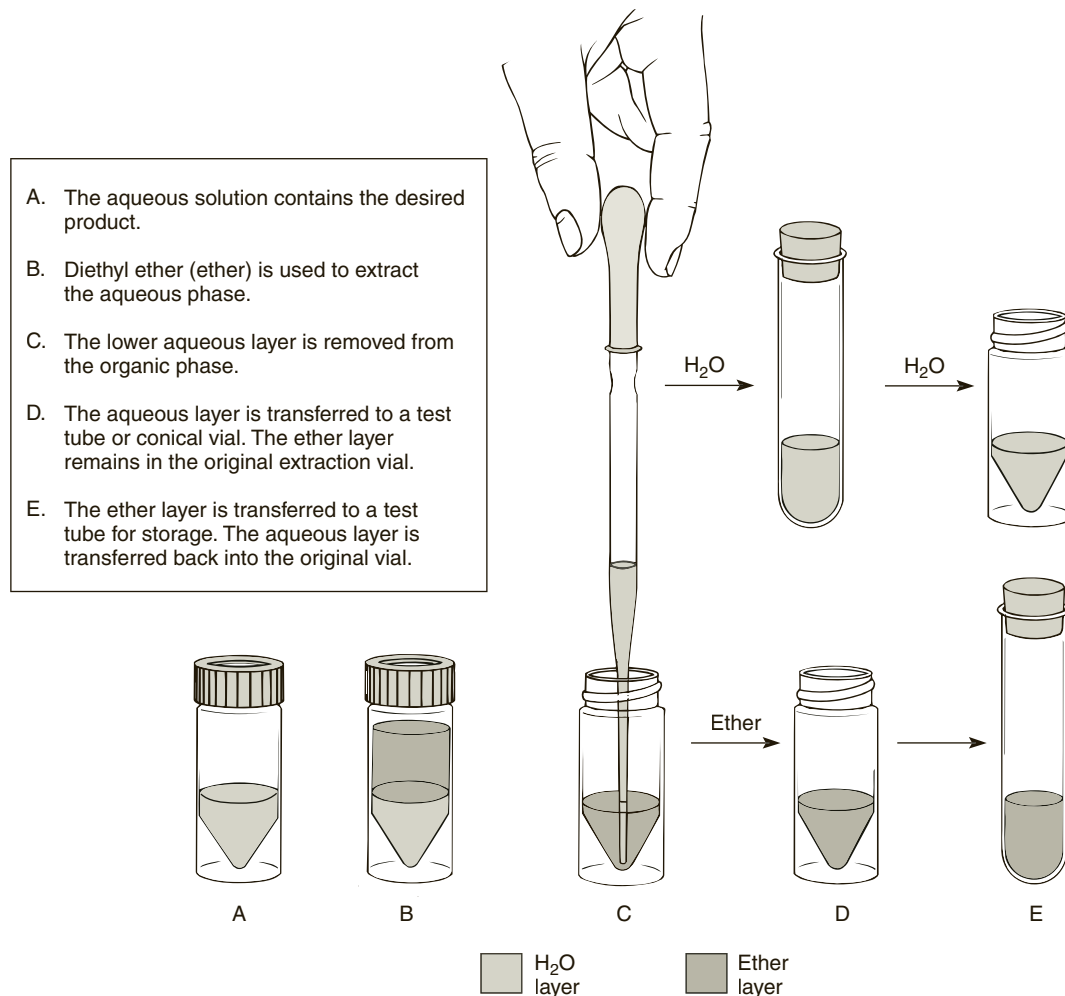**Figure 12.7**

Extraction of an aqueous solution using a solvent less dense than water: diethyl ether.

operation. With experience, you should be able to judge how much to squeeze the bulb to draw in the desired volume of liquid. Slowly draw the lower *aqueous* layer into the pipette. Be sure to keep the tip of the pipette squarely in the V at the bottom of the vial (Figure 12.7C).

4. Transfer the withdrawn aqueous phase into a test tube or another conical vial for temporary storage. It is best to have the test tube or vial located next to the extraction vial. This avoids messy and disastrous transfers. Hold the vials in the same hand between your index finger and thumb as shown in Figure 12.6. The ether layer is left behind in the conical vial (Figure 12.7D).
5. The ether phase remaining in the original conical vial should be transferred with a Pasteur pipette into a test tube for storage and the aqueous phase returned to the original conical vial (Figure 12.7E).

In performing an actual extraction, you would extract the aqueous phase with another 1-mL portion of fresh ether to achieve a more complete extraction. Steps 2–5 would be repeated, and the organic layers from both extractions would be combined in the test tube. In some cases, you may need to extract the aqueous layer a third time with yet another 1-mL portion of ether. Again, the ether would be combined with the other two layers. This overall process uses three 1-mL portions of

ether to transfer the product from the water layer into ether. The ether extracts contain some water and must be dried with a drying agent as indicated in Section 12.9.

12.6 The Centrifuge Tube

A screw-cap centrifuge tube may be employed instead of a conical vial for separations (Figure 12.4). Before using the centrifuge tube, be sure to check it for leaks as indicated in Section 12.3. You should use the same extraction and separation procedures described in Sections 12.4 and 12.5. You may also use a “regular” (nonscrew-cap) centrifuge tube for extractions, although it will be necessary to cork the tube before shaking it. Because a regular centrifuge tube will probably leak around the cork, it is best to mix the contents with a vortex mixer (Section 12.3) to avoid shaking the tube. If an emulsion has formed after mixing or shaking, you can use a centrifuge to aid in the separation of the layers (Section 12.10). Once the layers have separated, it is easy to use a Pasteur pipette to withdraw the lower layer from the tapered bottom of the centrifuge tube.

PART C. MACROSCALE EXTRACTION

12.7 The Separatory Funnel

The separatory funnel is often used in large-scale reactions. This apparatus is illustrated in Figure 12.8. To fill the separatory funnel, support it in an iron ring attached to a ring stand. Cut pieces of rubber tubing and attach them to the iron ring to cushion the separatory funnel as shown in Figure 12.8. This protects the funnel against possible breakage.

When beginning an extraction, the first step is to close the stopcock. (Don't forget!) Using a powder funnel (wide bore) placed in the top of the separatory funnel, fill it with both the solution to be extracted and the extraction solvent. Swirl the funnel gently by holding it by its upper neck, and then stopper it. Pick up the separatory funnel with two hands and hold it as shown in Figure 12.9. Hold the stopper in place firmly because the two immiscible liquids will build pressure when they mix, and this pressure may force the stopper out of the separatory funnel. To release this pressure, vent the funnel by holding it upside down (hold the stopper securely) and slowly open the stopcock. Usually the rush of vapors out of the opening can be heard. Continue shaking and venting until the “whoosh” is no longer audible. Now continue shaking the mixture gently for about one minute. This can be done by inverting the funnel in a rocking motion repeatedly or, if the formation of an emulsion is not a problem (see Section 12.10), by shaking the funnel more vigorously for less time.

NOTE: There is an art to shaking and venting a separatory funnel correctly, and it usually seems awkward to the beginner. The technique is best learned by observing a person, such as your instructor, who is thoroughly familiar with the separatory funnel's use.

When you have finished mixing the liquids, place the separatory funnel in the iron ring and remove the top stopper immediately. The two immiscible solvents separate into two layers after a short time, and they can be separated from each other by draining most of the lower layer through the stopcock.¹ Allow a few minutes to pass so that any of the lower phase adhering to the inner glass surfaces of the separatory funnel can drain

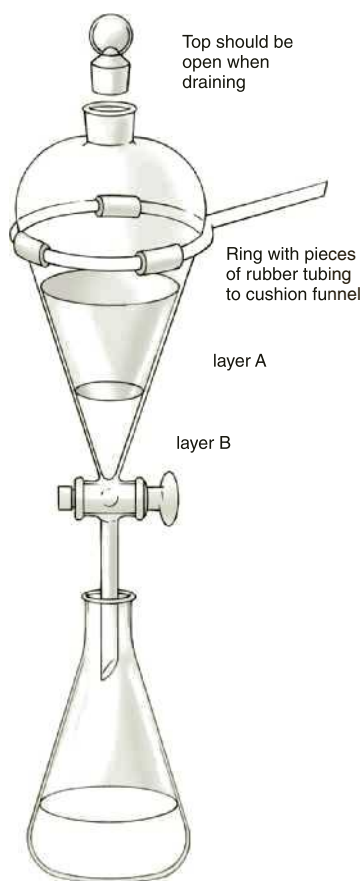


Figure 12.8
The separatory funnel.

¹ A common error is to try to drain the separatory funnel without removing the top stopper. Under this circumstance, the funnel will not drain, because a partial vacuum is created in the space above the liquid.

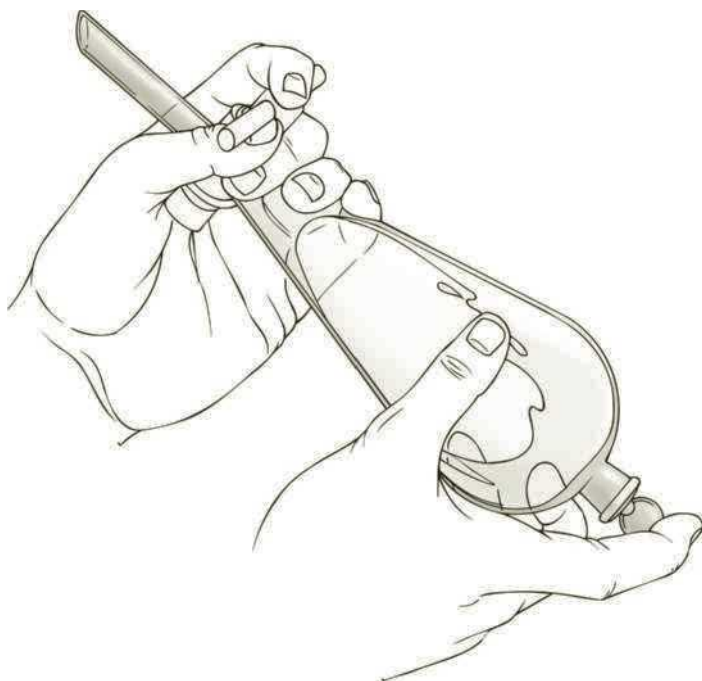


Figure 12.9
Correct way of shaking and venting the separatory funnel.

down. Open the stopcock again and allow the remainder of the lower layer to drain until the interface between the upper and lower phases just begins to enter the bore of the stopcock. At this moment, close the stopcock and remove the remaining upper layer by pouring it from the top opening of the separatory funnel.

NOTE: To minimize contamination of the two layers, the lower layer should always be drained from the bottom of the separatory funnel and the upper layer poured out from the top of the funnel.

When methylene chloride is used as the extracting solvent with an aqueous phase, it will settle to the bottom and be removed through the stopcock. The aqueous layer remains in the funnel. A second extraction of the remaining aqueous layer with fresh methylene chloride may be needed.

With a diethyl ether (ether) extraction of an aqueous phase, the organic layer will form on top. Remove the lower aqueous layer through the stopcock and pour the upper ether layer from the top of the separatory funnel. Pour the aqueous phase back into the separatory funnel and extract it a second time with fresh ether. The combined organic phases must be dried using a suitable drying agent (Section 12.9) before the solvent is removed.

For microscale procedures, a 60- or 125-mL separatory funnel is recommended. Because of surface tension, water has a difficult time draining from the bore of smaller funnels. Funnels larger than 125 mL are simply too large for microscale experiments, and a good deal of material is lost in “wetting” their surfaces.

PART D. ADDITIONAL EXPERIMENTAL CONSIDERATIONS: MICROSCALE AND MACROSCALE

12.8 How Do You Determine Which One Is the Organic Layer?

A common problem encountered during an extraction is trying to determine which of the two layers is the organic layer and which is the aqueous (water) layer. The most common situation occurs when the aqueous layer is on the bottom in the presence of an upper organic layer consisting of ether, ligroin, petroleum ether, or hexane (see densities in Table 12.1). However, the aqueous layer will be on the top when you use methylene chloride as a solvent (again, see Table 12.1). Although a laboratory procedure may frequently identify the expected relative positions of the organic and aqueous layers, sometimes their actual positions are reversed. Surprises usually occur in situations in which the aqueous layer contains a high concentration of sulfuric acid or a dissolved ionic compound, such as sodium chloride. Dissolved substances greatly increase the density of the aqueous layer, which may lead to the aqueous layer being found on the bottom even when coexisting with a relatively dense organic layer such as methylene chloride.


NOTE: Always keep both layers until you have actually isolated the desired compound or until you are certain where your desired substance is located.

To determine if a particular layer is the aqueous one, add a few drops of water to the layer. Observe closely as you add the water to see where it goes. If the layer is water, then the drops of added water will dissolve in the aqueous layer and increase its volume. If the added water forms droplets or a new layer, however, you can assume that the suspected aqueous layer is actually organic. You can use a similar procedure to identify a suspected organic layer. This time, try adding more of the solvent, such as methylene chloride. The organic layer should increase in size, without separation of a new layer, if the tested layer is actually organic.

When performing an extraction procedure on the microscale level, you can use the following approach to identify the layers. When both layers are present, it is always a good idea to think carefully about the volumes of materials that you have added to the conical vial. You can use the graduations on the vial to help determine the volumes of the layers in the vial. If, for example, you have 1 mL of methylene chloride in a vial and you add 2 mL of water, you should expect the water to be on top because it is less dense than methylene chloride. As you add the water, watch to see where it goes. By noting the relative volumes of the two layers, you should be able to tell which is the aqueous layer and which is the organic layer. This approach can also be used when performing an extraction procedure using a centrifuge tube. Of course, you can always test to see which layer is the aqueous layer by adding one or two drops of water, as described previously.

12.9 Drying Agents

After an organic solvent has been shaken with an aqueous solution, it will be “wet”; that is, it will have dissolved some water even though its solubility with water is not great. The amount of water dissolved varies from solvent to solvent; diethyl ether represents a solvent in which a fairly large amount of water dissolves. To remove water from the organic layer, use a **drying agent**. A drying agent is an *anhydrous* inorganic salt that acquires waters of hydration when exposed to moist air or a wet solution:

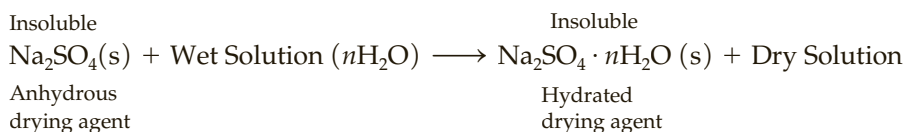

The insoluble drying agent is placed directly into the solution, where it acquires water molecules and becomes hydrated. If enough drying agent is used, all of the water can be removed from a wet solution, making it “dry,” or free of water.

The following anhydrous salts are commonly used: sodium sulfate, magnesium sulfate, calcium chloride, calcium sulfate (Drierite), and potassium carbonate. These salts vary in their properties and applications. For instance, not all will absorb the same amount of water for a given weight, nor will they dry the solution to the same extent. **Capacity** refers to the amount of water a drying agent absorbs per unit weight. Sodium and magnesium sulfates absorb a large amount of water (high capacity), but magnesium sulfate dries a solution more completely. **Completeness** refers to a compound’s effectiveness in removing all the water from a solution by the time equilibrium has been reached. Magnesium ion, a strong Lewis acid, sometimes causes rearrangements of compounds such as epoxides. Calcium chloride is a good drying agent but cannot be used with many compounds containing oxygen or nitrogen because it forms complexes. Calcium chloride absorbs methanol and ethanol in addition to water, so it is useful for removing these materials when they are present as impurities. Potassium carbonate is a base and is used for drying solutions of basic substances, such as amines. Calcium sulfate dries a solution completely but has a low capacity.

Table 12.2 Common drying agents

	<i>Acidity</i>	<i>Hydrated</i>	<i>Capacity</i> ^a	<i>Completeness</i> ^b	<i>Rate</i> ^c	<i>Use</i>
Magnesium sulfate	Neutral	MgSO ₄ ·7H ₂ O	High	Medium	Rapid	General
Sodium sulfate	Neutral	Na ₂ SO ₄ ·7H ₂ O Na ₂ SO ₄ ·10H ₂ O	High	Low	Medium	General
Calcium chloride	Neutral	CaCl ₂ ·2H ₂ O CaCl ₂ ·6H ₂ O	Low	High	Rapid	Hydrocarbons Halides
Calcium sulfate (Drierite)	Neutral	CaSO ₄ ·½H ₂ O CaSO ₄ ·2H ₂ O	Low	High	Rapid	General
Potassium carbonate	Basic	K ₂ CO ₃ ·1½H ₂ O K ₂ CO ₃ ·2H ₂ O	Medium	Medium	Medium	Amines, esters, bases, ketones
Potassium hydroxide	Basic	—	—	—	Rapid	Amines only
Molecular sieves (3 or 4 Å)	Neutral	—	High	Extremely high	—	General

^aAmount of water removed per given weight of drying agent.

^bRefers to amount of H₂O still in solution at equilibrium with drying agent.

^cRefers to rate of action (drying).

Anhydrous sodium sulfate is the most widely used drying agent. The granular variety is recommended because it is easier to remove the dried solution from it than from the powdered variety. Sodium sulfate is mild and effective. It will remove water from most common solvents, with the possible exception of diethyl ether, in which case a prior drying with saturated salt solution may be advised. Sodium sulfate must be used at room temperature to be effective; it cannot be used with boiling solutions. Table 12.2 compares the various common drying agents.

Drying Procedure with Anhydrous Sodium Sulfate. In experiments that require a drying step, the instructions are usually given in the following way: dry the organic layer (or phase) over granular anhydrous sodium sulfate (or some other drying agent). More specific instructions, such as the amount of drying agent to add, usually will not be given, and you will need to determine this each time that you perform a drying step. The drying procedure consists of four steps:

1. Remove the organic layer from any visible water.
2. Add the appropriate amount of granular anhydrous sodium sulfate (or other drying agent).
3. Allow a drying period during which dissolved water is removed from the organic layer by the drying agent.
4. Separate the dried organic layer from the drying agent.

More specific instructions are given below for both microscale and macroscale procedures. The only differences between these two procedures is that they are intended for different volumes of liquid and they require different glassware. The microscale procedure is generally for volumes up to about 5 mL, and the macroscale procedure is usually appropriate for volumes of 5 mL or greater.

A. Microscale Drying Procedure

Step 1. (Removal of Visible Water). Before attempting to dry an organic layer, check closely to see that there are no visible signs of water. If there is a separate layer of water (top or bottom), droplets or a globule of water floating in the organic layer, or water droplets clinging to the sides of the container, then transfer the organic layer with a *dry* Pasteur pipette to a *dry* container, usually a conical vial or test tube, before adding any drying agent. If there is any doubt about whether water is present, it is advisable to make a transfer to a dry container. Performing this step when necessary may save time later in the drying procedure and result in a greater recovery of the desired substance.

Step 2. (Addition of Drying Agent). Each time a drying procedure is performed, it is necessary to determine how much granular anhydrous sodium sulfate (or other drying agent) should be added. This will depend on the total volume of the organic phase and how much water is dissolved in the solvent. Nonpolar organic solvents such as methylene chloride or hydrocarbons (hexane, pentane, etc.) can dissolve relatively small amounts of water and generally require less drying agent, whereas more polar organic solvents such as ether and ethyl acetate can dissolve more water, and more drying agent will be required. Begin by adding one spatulaful of granular anhydrous sodium sulfate (or other drying agent) from the V-grooved end of a microspatula (smaller microspatula in Experiment 1, Figure 10) into the solution. If all the drying agent “clumps,” add another spatulaful of sodium sulfate. To determine if the drying agent has clumped, it is helpful to stir the mixture with a clean, dry spatula or to rapidly swirl the container. If any portion of the drying agent flows freely (does not clump) on the bottom of the container when stirred or swirled, then you can assume that enough of the drying agent has been added. Otherwise, you must continue adding one spatulaful of drying agent at a time until it is clear that the drying agent has stopped clumping. Stir or swirl the mixture after adding each spatulaful of the drying agent. For small amounts of liquid (less than 5 mL), about 1–6 microspatulafuls of drying agent will *usually* be required. However, the actual amount must be determined by experimentation, as just described. It is best to use a slight excess of drying agent; but if too great an excess is used, the recovery may be poor because some of the solution always adheres to the solid drying agent after the liquid is separated from the drying agent (Step 4). Take care not to add so much drying agent that all of the liquid is absorbed (disappears). If you do this you will have to add additional solvent to recover your product from the drying agent!

Step 3. (Drying Period). Stopper or cap the container and let the solution dry for at least 15 minutes.

NOTE: It is important that you stopper or cap the container to prevent evaporation and exposure to atmospheric moisture.

Stir the mixture occasionally with a spatula during the drying period. The mixture is dry if it appears clear (not cloudy) and shows the common signs of a dry solution given in Table 12.3. Note that a “clear” solution may be colorless or colored. If the solution remains cloudy after treatment with the first batch of drying agent, add more drying agent and repeat the drying procedure. However, if a water layer forms or if drops of water are visible, transfer the organic layer to a dry container before adding fresh drying agent, as described in Step 2. It will also be necessary to repeat the 15-minute drying step described in Step 3.

Table 12.3 Common signs that indicate a solution is dry

1. There are no visible water droplets on the side of flask or suspended in solution.
2. There is not a separate layer of liquid or a “puddle.”
3. The solution is clear, not cloudy. Cloudiness indicates water is present.
4. The drying agent (or a portion of it) flows freely on the bottom of the container when stirred or swirled and does not “clump” together as a solid mass.

Step 4. (Removal of Liquid from Drying Agent). When the organic phase is dry, use a *dry* Pasteur pipette or a *dry* filter-tip pipette (Technique 8, Section 8.6) to remove the dried organic layer from the drying agent and transfer the solution to a *dry* conical vial or test tube. Be careful not to transfer any of the drying agent when performing this step. Rinse the drying agent with a small amount of fresh solvent and transfer this additional solvent to the vial containing the dried organic layer. To isolate the desired material, remove the solvent by evaporation using heat and a stream of air or nitrogen (Technique 7, Section 7.10).

An alternative method of drying a small volume of organic phase is to pass it through a filtering pipette (Technique 8, Section 8.1C) that has been packed with a small amount (about 2 cm) of drying agent. Again, the solvent is removed by evaporation.

B. Macroscale Drying Procedure

To dry a large amount of organic liquid (greater than about 5 mL), follow the same four steps just described for the “Microscale Drying Procedure.” The main differences are that an Erlenmeyer flask is used rather than a test tube or conical vial and more drying agent will be required. The size of the Erlenmeyer flask is not critical, but it’s best that the flask not be filled more than half full with the solution being dried.

Step 1. (Removal of Visible Water). Refer to Step 1 above for instructions. If the amount of water is large, it may be best to separate the layers using a separatory funnel or a centrifuge tube. If visible water must be removed in this step, place the separated organic layer in a clean, *dry* Erlenmeyer flask.

Step 2. (Addition of Drying Agent). Refer to Step 2 in the “Microscale Drying Procedure” for the basic instructions. Read these instructions carefully. The only difference is that in this macroscale procedure, more drying agent will be required. A common guideline is to add enough granular anhydrous sodium sulfate (or other drying agent) to give a 1- to 3-mm layer on the bottom of the flask, depending on the volume of the solution. However, it is best to add the drying agent in small portions, as described above. In this procedure, use the larger microspatula shown in Experiment 1, Figure 10 to add the drying agent. Generally, an appropriate portion to add each time is about 0.5–1.0 g.

Step 3. (Drying Period). The instructions are the same as for Step 3 in the “Microscale Drying Procedure.”

Step 4. (Removal of Liquid from Drying Agent). When the solution is dry, the drying agent should be removed by using decantation (pouring carefully to leave the drying agent behind). Transfer the liquid to a *dry* Erlenmeyer flask. If the volume of liquid is relatively small (less than 10 mL), it may be easier to complete this step by using a *dry* Pasteur pipette or a *dry* filter-tip pipette (Technique 8, Section 8.6)

to remove the dried organic layer. With granular sodium sulfate, decantation is easy to perform because of the size of the drying-agent particles. If a powdered drying agent, such as magnesium sulfate, is used, it may be necessary to use gravity filtration (Technique 8, Section 8.1B) to remove the drying agent. Finally, to isolate the desired material, remove the solvent by distillation (Technique 14, Section 14.3) or evaporation (Technique 7, Section 7.10).

Saturated Salt Solution. At room temperature, diethyl ether (ether) dissolves 1.5% by weight of water, and water dissolves 7.5% of ether. Ether, however, dissolves a much smaller amount of water from a saturated aqueous sodium chloride solution. Hence, the bulk of water in ether, or ether in water, can be removed by shaking it with a saturated aqueous sodium chloride solution. A solution of high ionic strength is usually not compatible with an organic solvent and forces separation of it from the aqueous layer. The water migrates into the concentrated salt solution. The ether phase (organic layer) will be on top, and the saturated sodium chloride solution will be on the bottom ($d = 1.2 \text{ g/mL}$). After removing the organic phase from the aqueous sodium chloride, dry the organic layer completely with sodium sulfate or with one of the other drying agents listed in Table 12.2.

12.10 Emulsions

An **emulsion** is a colloidal suspension of one liquid in another. Minute droplets of an organic solvent are often held in suspension in an aqueous solution when the two are mixed or shaken vigorously; these droplets form an emulsion. This is especially true if any gummy or viscous material was present in the solution. Emulsions are often encountered in performing extractions. Emulsions may require a long time to separate into two layers and are a nuisance to the organic chemist.

Fortunately, several techniques may be used to break a difficult emulsion once it has formed.

1. Often an emulsion will break up if it is allowed to stand for some time. Patience is important here. Gently stirring with a stirring rod or spatula may also be useful.
2. If one of the solvents is water, adding a saturated aqueous sodium chloride solution will help destroy the emulsion. The water in the organic layer migrates into the concentrated salt solution.
3. With microscale experiments, the mixture may be transferred to a centrifuge tube. The emulsion will often break during centrifugation. Remember to place another tube filled with water on the opposite side of the centrifuge to balance it. The two tubes should weigh the same.
4. Adding a small amount of a water-soluble detergent may also help. This method has been used in the past for combating oil spills. The detergent helps to solubilize the tightly bound oil droplets.
5. Gravity filtration (see Technique 8, Section 8.1) may help to destroy an emulsion by removing gummy polymeric substances. With large volumes, you might try filtering the mixture through a fluted filter (Technique 8, Section 8.1B) or a piece of cotton. With small-scale reactions, a filtering pipette may work (Technique 8, Section 8.1C). In many cases, once the gum is removed, the emulsion breaks up rapidly.
6. If you are using a separatory funnel, you might try to use a gentle swirling action in the funnel to help break an emulsion. Gently stirring with a stirring rod may also be useful.


When you know through experience that a mixture may form a difficult emulsion, you should avoid shaking the mixture vigorously. When using conical vials for extractions, it may be better to use a magnetic spin vane for mixing and not shake the mixture at all. When using separatory funnels, extractions should be performed with gentle swirling instead of shaking or with several gentle inversions of the separatory funnel. Do not shake the separatory funnel vigorously in these cases. It is important to use a longer extraction period if the more gentle techniques described in this paragraph are being employed. Otherwise, you will not transfer all the material from the first phase to the second one.

12.11 Purification and Separation Methods

In nearly all synthetic experiments undertaken in the organic laboratory, a series of operations involving extractions is used after the actual reaction has been concluded. These extractions form an important part of the purification. Using them, you separate the desired product from unreacted starting materials or from undesired side products in the reaction mixture. These extractions may be grouped into three categories, depending on the nature of the impurities they are designed to remove.


The first category involves extracting or “washing” an organic mixture with water. Water washes are designed to remove highly polar materials, such as inorganic salts, strong acids or bases, and low-molecular-weight, polar substances including alcohols, carboxylic acids, and amines. Many organic compounds containing fewer than five carbons are water soluble. Water extractions are also used immediately following extractions of a mixture with either acid or base to ensure that all traces of acid or base have been removed.

The second category concerns extraction of an organic mixture with a dilute acid, usually 1–2 M hydrochloric acid. Acid extractions are intended to remove basic impurities, especially such basic impurities as organic amines. The bases are converted to their corresponding cationic salts by the acid used in the extraction. If an amine is one of the reactants or if pyridine or another amine is a solvent, such an extraction might be used to remove any excess amine present at the end of a reaction.

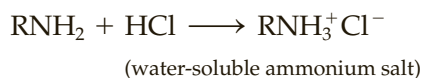


Cationic ammonium salts are usually soluble in the aqueous solution, and they are thus extracted from the organic material. A water extraction may be used immediately following the acid extraction to ensure that all traces of the acid have been removed from the organic material.

The third category is extraction of an organic mixture with a dilute base, usually 1 M sodium bicarbonate, although extractions with dilute sodium hydroxide can also be used. Such basic extractions are intended to convert acidic impurities, such as organic acids, to their corresponding anionic salts. For example, in the preparation of an ester, a sodium bicarbonate extraction might be used to remove any excess carboxylic acid that is present.

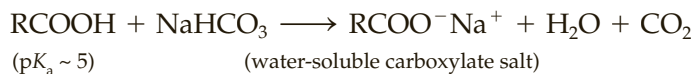

Anionic carboxylate salts, being highly polar, are soluble in the aqueous phase. As a result, these acid impurities are extracted from the organic material into the basic solution. A water extraction may be used after the basic extraction to ensure that all the base has been removed from the organic material.

Occasionally, phenols may be present in a reaction mixture as impurities, and removing them by extraction may be desired. Because phenols, although they are acidic, are about 10^5 times less acidic than carboxylic acids, basic extractions may be used to separate phenols from carboxylic acids by a careful selection of the base. If sodium bicarbonate is used as a base, carboxylic acids are extracted into the aqueous base, but phenols are not. Phenols are not sufficiently acidic to be deprotonated by the weak base bicarbonate. Extraction with sodium hydroxide, on the other hand, extracts both carboxylic acids and phenols into the aqueous basic solution because hydroxide ion is a sufficiently strong base to deprotonate phenols.

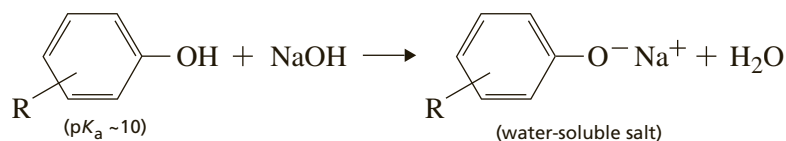

Mixtures of acidic, basic, and neutral compounds are easily separated by extraction techniques. One such example is shown in Figure 12.10. The original compounds are dissolved in ether.

Organic acids or bases that have been extracted can be regenerated by neutralizing the extraction reagent. This would be done if the organic acid or base were a product of a reaction rather than an impurity. For example, if a carboxylic acid has been extracted with the aqueous base, the compound can be regenerated by acidifying the extract with 6 M HCl until the solution becomes *just* acidic, as indicated by litmus or pH paper. When the solution becomes acidic, the carboxylic acid will separate from the aqueous solution. If the acid is a solid at room temperature, it will precipitate and can be purified by filtration and crystallization. If the acid is a liquid, it will form a separate layer. In this case, it would usually be necessary to

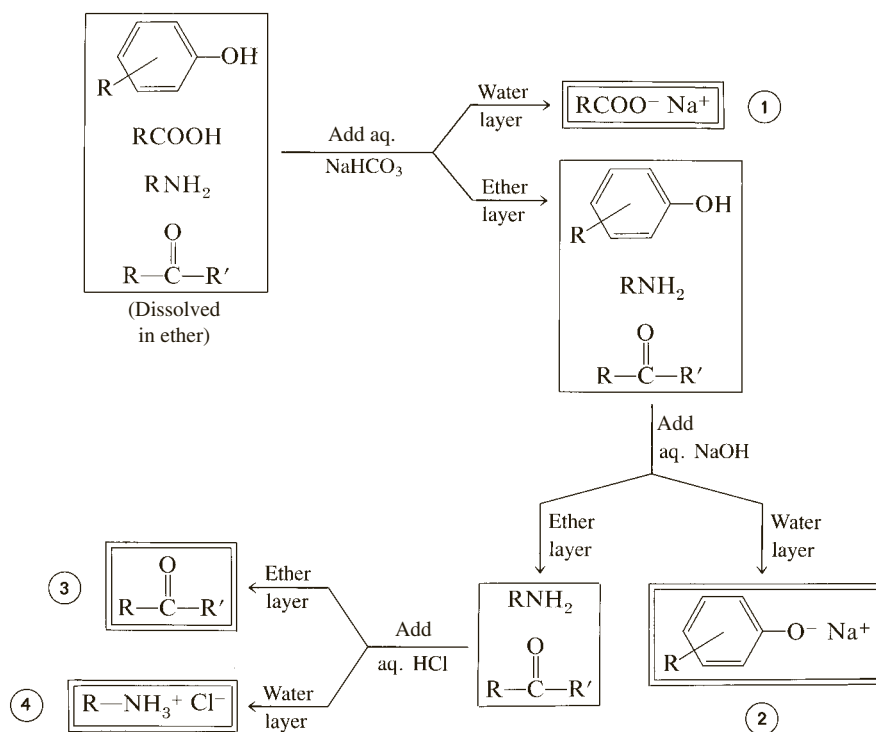


Figure 12.10
Separating a four-component mixture by extraction.

extract the mixture with ether or methylene chloride. After removing the organic layer and drying it, the solvent can be evaporated to yield the carboxylic acid.

In the example shown in Figure 12.10, you also need to perform a drying step at (3) before isolating the neutral compound. When the solvent is ether, you should first extract the ether solution with saturated aqueous sodium chloride to remove much of the water. The ether layer is then dried over a drying agent such as anhydrous sodium sulfate. If the solvent were methylene chloride, it would not be necessary to do the step with saturated sodium chloride.

When acid–base extractions are performed, it is common practice to extract a mixture several times with the appropriate reagent. For example, if you were extracting a carboxylic acid from a mixture, you might extract the mixture three times with 2-mL portions of 1 M NaOH. In most published experiments, the procedure will specify the volume and concentration of extracting reagent and the number of times to do the extractions. If this information is not given, you must devise your own procedure. Using a carboxylic acid as an example, if you know the identity of the acid and the approximate amount present, you can actually calculate how much sodium hydroxide is needed. Because the carboxylic acid (assuming it is monoprotic) will react with sodium hydroxide in a 1:1 ratio, you would need the same number of moles of sodium hydroxide as there are moles of acid. To ensure that all the carboxylic acid is extracted, you should use about a *twofold* excess of the base. From this, you could calculate the number of milliliters of base needed. This should be divided into two or three equal portions, one portion for each extraction. In a similar fashion, you could calculate the amount of 5% sodium bicarbonate required to extract an acid or the amount of 1 M HCl required to extract a base. If the amount of organic acid or base is not known, then the situation is more difficult. A guideline that sometimes works is to do two or three extractions so that the total volume of the extracting reagent is approximately equal to the volume of the organic layer. To test this procedure, neutralize the aqueous layer from the last extraction. If a precipitate or cloudiness results, perform another extraction and test again. When no precipitate forms, you know that all the organic acid or base has been removed.

For some applications of acid–base extraction, an additional step, called **backwashing** or **back extraction**, is added to the scheme shown in Figure 12.10. Consider the first step, in which the carboxylic acid is extracted by sodium bicarbonate. This aqueous layer may contain some unwanted neutral organic material from the original mixture. To remove this contamination, backwash the aqueous layer with an organic solvent such as ether or methylene chloride. After shaking the mixture and allowing the layers to separate, remove and discard the organic layer. This technique may also be used when an amine is extracted with hydrochloric acid. The resulting aqueous layer is backwashed with an organic solvent to remove unwanted neutral material.

12.12 How to Outline a Separation Scheme Using a Flowchart

In experiments involving the synthesis of an organic compound, it is necessary to separate the desired product from the reaction mixture, which usually contains many other substances. The steps required to accomplish this are sometimes referred to as a separation scheme. An outline of a separation scheme can be conveniently shown in the form of a flowchart. A flowchart for the separation scheme for isopentyl acetate is shown in Technique 2, Figure 2.1. Experimental techniques that may be included in a separation scheme include acid–base extraction, a drying step, sublimation, distillation, and crystallization. It is particularly helpful to outline a separation scheme before coming to the laboratory. Preparing the flowchart will help you to understand all of the steps in the experiment, which may lead to greater success when actually performing an experiment.

The starting point of most separation schemes is at the end of the reaction period. Do not include any of the reaction steps such as adding reagents or the reflux period. The reaction mixture will often contain the desired product, leftover reactants, catalysts, solvents, by-products that are also formed in the main reaction, and side products that are produced by undesired reactions that take place along with the main reaction. The purpose of a separation scheme is to separate the desired product from all other substances, resulting in a pure product. At the beginning of your flowchart, you should list everything that is present at the end of the reaction period, either by name or formula. It may sometimes be convenient to use a condensed or partial formula. For example, octanoic acid might be written as $C_8H_{17}COOH$ or just $R-COOH$, if this is defined in the beginning list. Every step involved in separating the desired product from the other substances should be shown in your flowchart. When a new reagent such as hydrochloric acid, sodium hydroxide, or anhydrous sodium sulfate is added, you should indicate this in your flowchart and then you must also show how this reagent is removed in the separation scheme. You should indicate by name any experimental procedure that is used, such as distillation, sublimation, and crystallization. If you don't know the physical properties of the product, you should look this up in a handbook before writing the flowchart. This is important because a solid product would require different techniques to separate than a liquid product.


It would be helpful now to study the separation scheme for isopentyl acetate shown in Technique 2, Figure 2.1. Note carefully how the following procedures are shown in a flowchart: the extraction step with $NaHCO_3$, the drying step with anhydrous Na_2SO_4 , and the distillation step. Also note that in this experiment it is not possible to identify with certainty the actual impurities present with isopentyl acetate just before the product is distilled. In some experiments it will be possible to make a good guess about the identity of most substances at each step in the separation scheme.

There are several other experimental techniques that may come up in separation schemes that are not shown in Figure 2.1. We will now discuss some of these additional techniques and show how to use a flowchart to illustrate these techniques.


First, consider the extraction of a solid carboxylic acid (the desired product) with sodium bicarbonate, followed by precipitation of the acid and isolation of the solid acid by filtration. It is also necessary to air-dry the filtered solid. See the flowchart for all of these steps in Figure 12.11. Although this procedure sometimes gives the acid in a pure form, it may also be necessary to perform a crystallization procedure. In this example, the carboxylic acid is assumed to be a solid at room temperature. If the carboxylic acid were a liquid at room temperature, the separation scheme would involve different techniques including a distillation. The neutral impurities are not identified here since this example represents a general case rather than a specific experiment. In many experiments it will be possible to determine the identity of most of the impurities. For example, in most synthesis experiments, some of the impurities will be leftover reactants.

Note that the use of $NaHCO_3$ in this separation scheme is different from the $NaHCO_3$ extraction step in Figure 2.1. In the separation scheme shown in Figure 2.1, the $NaHCO_3$ step is used to remove acetic acid from the mixture containing the product and it is not necessary to regenerate the acetic acid in a subsequent step. However, in the procedure shown in Figure 12.11, the carboxylic acid is the desired product and it must be regenerated as a precipitate by the addition of HCl so that it can be isolated.

In the second example, an amine, which is a solid base and the desired product, is separated from the reaction mixture in a similar set of steps, except that HCl is used to extract the base and $NaOH$ is used to regenerate the base. See Figure 12.12 for a flowchart of this separation scheme. The amine is represented by $R-NH_2$.

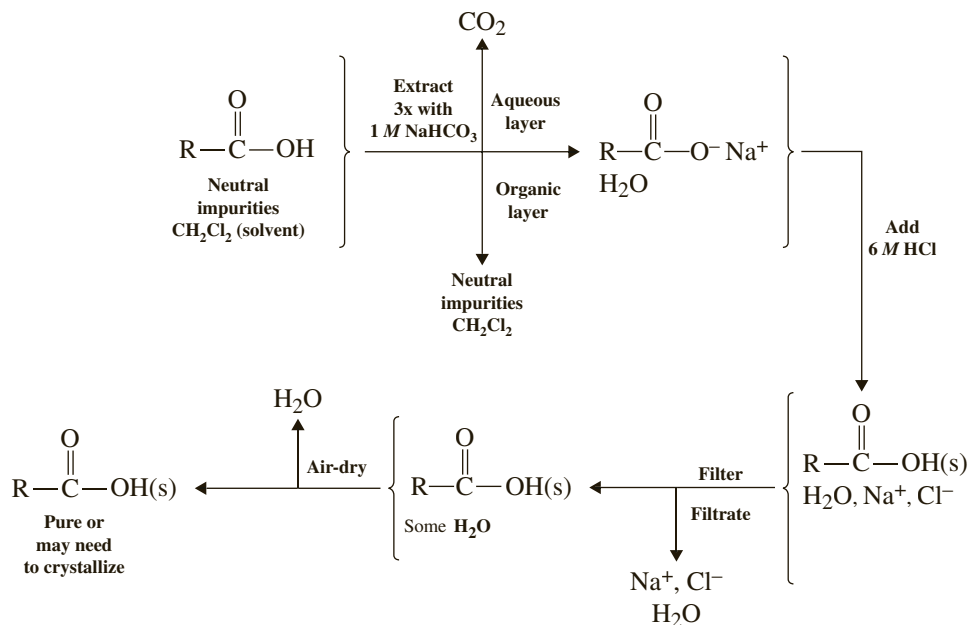**Figure 12.11**

Separation of a solid carboxylic acid from a reaction mixture.

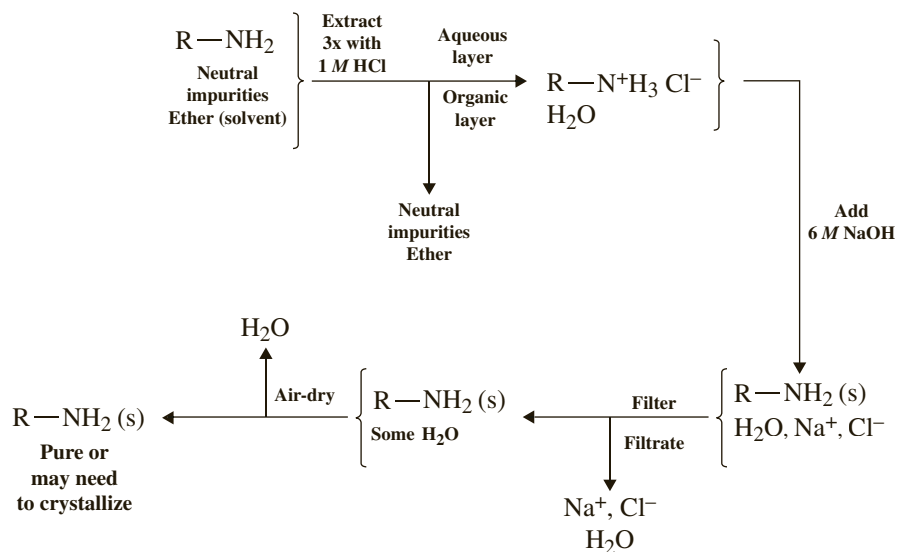**Figure 12.12**

Separation of a solid amine from a reaction mixture.

The last example involves a crystallization step in which a solid neutral compound is purified from an impure sample consisting mainly of the neutral compound and a second neutral compound which is an impurity found in a much smaller amount than the compound of interest. For example, in Experiment 6C, fluorenone is synthesized from fluorenone. The crude product contains mainly fluorenone contaminated with a little fluorenone. See Figure 12.13 for a flowchart showing this crystallization step.

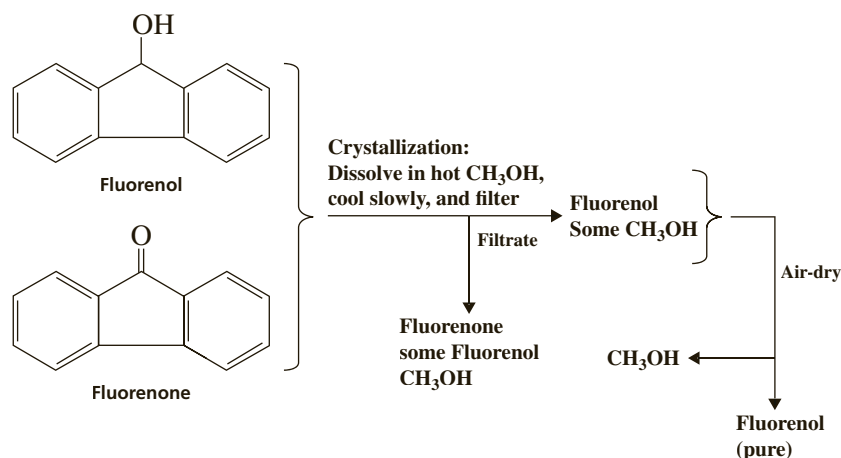


Figure 12.13
Purification of a neutral compound by crystallization.

PART E. OTHER EXTRACTION METHODS

12.13 Continuous Solid–Liquid Extraction

The technique of liquid–liquid extraction was described in Sections 12.1–12.8. In this section, solid–liquid extraction is described. Solid–liquid extraction is often used to extract a solid natural product from a natural source, such as a plant. A solvent is chosen that selectively dissolves the desired compound but that leaves behind the undesired insoluble solid. A continuous solid–liquid extraction apparatus, called a Soxhlet extractor, is commonly used in a research laboratory.

As shown in Figure 12.14, the solid to be extracted is placed in a thimble made from filter paper, and the thimble is inserted into the central chamber. A low-boiling solvent, such as diethyl ether, is placed in the round-bottom distilling flask and is heated to reflux. The vapor rises through the left side arm into the condenser where it liquefies. The condensate (liquid) drips into the thimble containing the solid. The hot solvent begins to fill the thimble and extracts the desired compound from the solid. Once the thimble is filled with solvent, the side arm on the right acts as a siphon, and the solvent, which now contains the dissolved compound, drains back into the distillation flask. The vaporization–condensation–extraction–siphoning process is repeated hundreds of times, and the desired product is concentrated in the distillation flask. The product is concentrated in the flask because the product has a boiling point higher than that of the solvent or because it is a solid.

12.14 Continuous Liquid–Liquid Extraction

When a product is very soluble in water, it is often difficult to extract using the techniques described in Sections 12.4–12.7 because of an unfavorable distribution coefficient. In this case, you need to extract the aqueous solution numerous times with fresh batches of an immiscible organic solvent to remove the desired product from water. A less labor-intensive technique involves the use of a continuous liquid–liquid extraction apparatus. One type of extractor, used with solvents that are less dense than water, is shown in Figure 12.15. Diethyl ether is usually the solvent of choice.

The aqueous phase is placed in the extractor, which is then filled with diethyl ether up to the side arm. The round-bottom distillation flask is partially filled with ether. The ether is heated to reflux in the round-bottom flask, and the vapor is liquefied in the water-cooled condenser. The ether drips into the central tube, passes

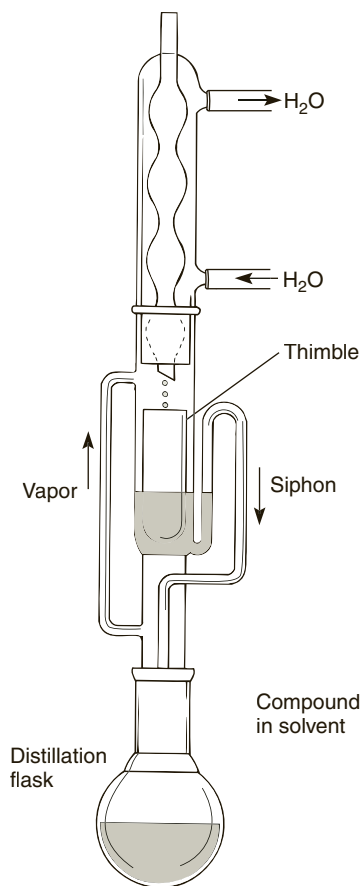


Figure 12.14
Continuous solid–liquid
extraction using a Soxhlet
extractor.

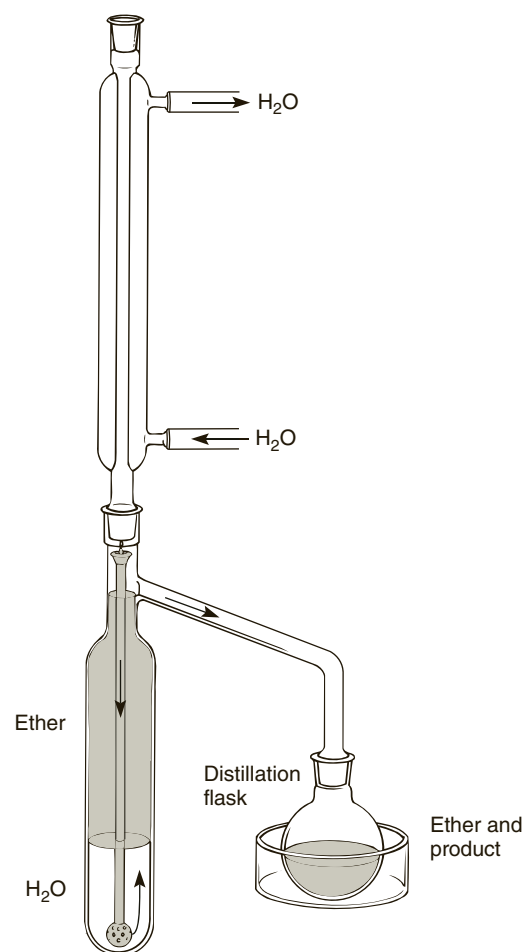


Figure 12.15
Continuous liquid–liquid extraction using a
solvent less dense than water.

through the porous sintered glass tip, and flows through the aqueous layer. The solvent extracts the desired compound from the aqueous phase, and the ether is recycled back into the round-bottom flask. The product is concentrated in the flask. The extraction is rather inefficient and must be placed in operation for at least 24 hours to remove the compound from the aqueous phase.

12.15 Solid Phase Extraction

Solid phase extraction (SPE) is a relatively new technique, which is similar in appearance and function to column chromatography and high-performance liquid chromatography (Techniques 19 and 21). In some applications, SPE is also similar to liquid–liquid extraction, discussed in this technique chapter. In addition to performing separation processes, SPE can also be used to carry out reactions in which new compounds are prepared.

A typical SPE column is constructed from the body of a plastic syringe, which is packed with a **sorbent**. The term *sorbent* is used by many manufacturers as a general term for materials that can both adsorb (attract to the surface of the sorbent by a physical attraction) or absorb (penetrate into the material like a sponge). A frit is inserted at the bottom of the column to support the sorbent. After the sorbent is added, another frit is inserted on top of the sorbent to hold it in place. The remainder of the tube serves as a reservoir for the solvent. Generally, the column comes packed with

the sorbent from the manufacturer, but unpacked columns can also be purchased and packed by the user for specific applications. The Luer-lock tip at the bottom is connected to a vacuum source that pulls the solvents through the column.

SPE columns can be packed with many kinds of sorbent, depending on how the column will be used. Some common types are identified in the same way that column chromatography adsorbents are classified (see Technique 21, Section 21.1): normal-phase, reversed-phase, and ion exchange. Examples of normal-phase sorbents, which are polar, include silica and alumina. These columns are used to isolate polar compounds from a nonpolar solvent. Reversed-phase sorbents are made by alkylating silica. As a result, nonpolar alkyl groups are bonded to the silica surface, making the sorbent nonpolar. A common column of this type, known as a C_{18} column, is prepared by attaching an octadecyl ($-C_8H_{18}$) group to the silica surface (see Figure 12.16). C_{18} columns most likely function by an adsorption process.

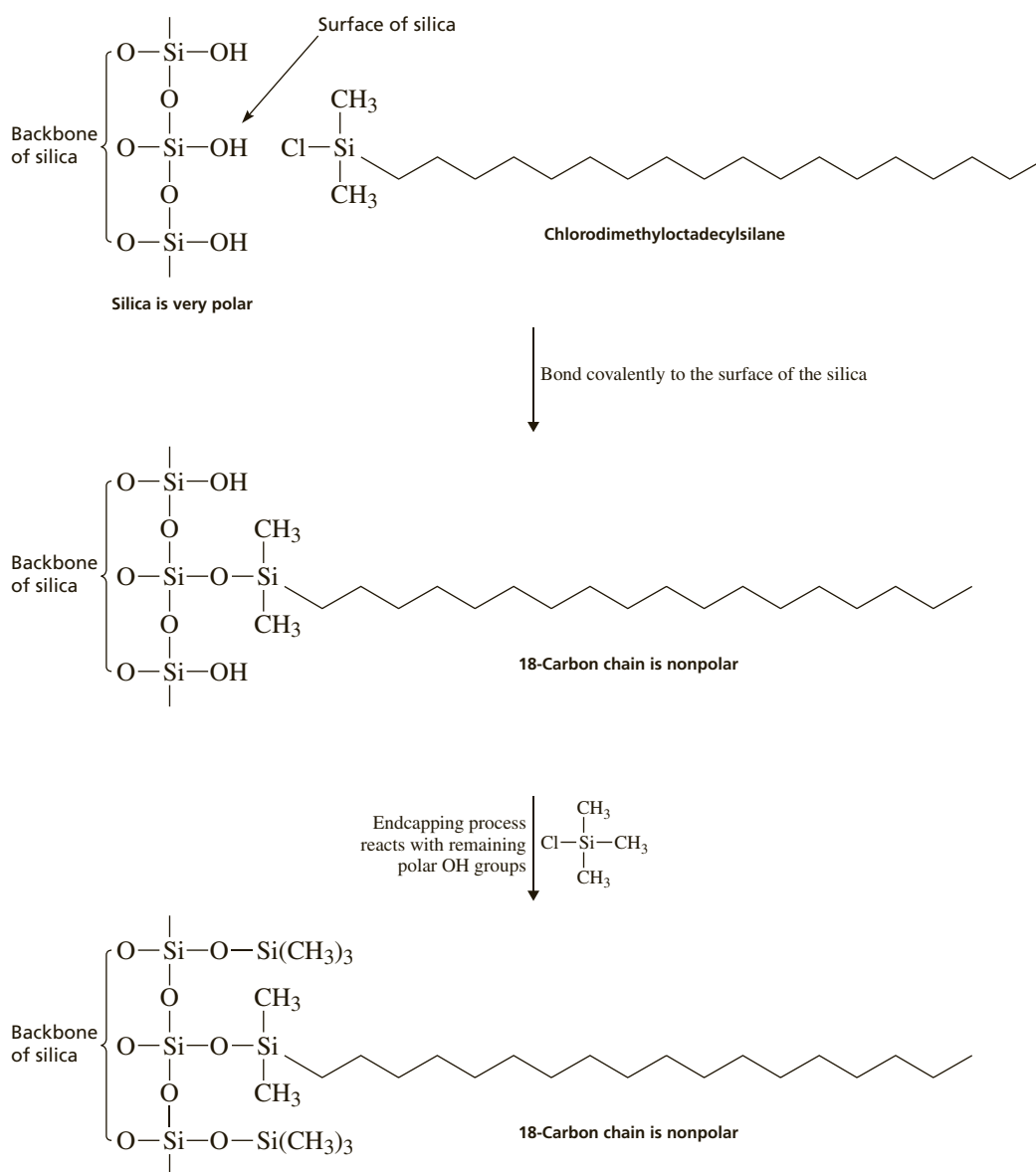


Figure 12.16

Preparation of C_{18} silica for reversed-phase extractions using SPE tubes. The process changes polar silica (hydrophilic material) to nonpolar silica (hydrophobic material).

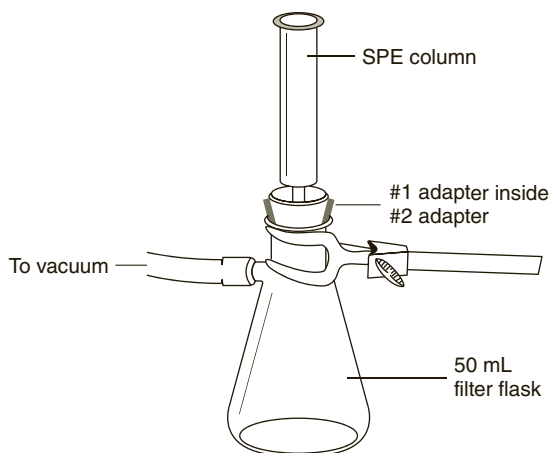


Figure 12.17
Experimental setup for SPE column.

Reversed-phase sorbents are used to isolate relatively nonpolar compounds from polar solvents. Ion-exchange sorbents consist of charged or highly polar materials and are used to isolate charged compounds, either as anions or cations.

A major advantage of SPE columns is that they are fast and convenient to use compared to traditional column chromatography or liquid–liquid extraction. However, there are many other advantages that are of benefit to the environment, and their use is a good example of green chemistry (see essay “Green Chemistry” that precedes Experiment 29). These advantages include the use of more environmentally friendly solvents, higher recovery, elimination of emulsions, enormous decrease in the use of solvents, and reduced toxic waste generation.

A good example of the use of SPE columns for performing a task that is normally done by liquid–liquid extraction is the isolation of caffeine from tea or coffee. In this application, a C_{18} column is used. As the tea or coffee flows through the column, caffeine is attracted to the sorbent, and the polar impurities come off with water. Ethyl acetate is then used to remove the caffeine from the column. The experimental setup for this is shown in Figure 12.17. The SPE column² is attached to the filter flask by using two neoprene adapters (sizes #1 and #2). The filter flask is connected to either a vacuum line or a water aspirator to provide the vacuum. After each step, the solvents with impurities or desired product are drawn through the column into the filter flask using the vacuum.

The following steps are used with an SPE tube to remove caffeine from tea or coffee (see Figure 12.18):

- A. Condition the C_{18} reversed-phase silica column by passing methanol and water through the tube.
- B. Apply the sample of caffeinated drink to the column.
- C. Wash the polar impurities from the column with water.
- D. Elute the caffeine from the tube with ethyl acetate.

Even though Figure 12.18 is applied to the isolation of caffeine, the general scheme may be used in any application in which it is desired to separate polar substances, such as water, from a relatively nonpolar substance. Numerous applications are found in the medical field, in which analyzing body fluids is important.

There are many other diverse applications that SPE columns can be used for. By modifying the silica with specific chemical reagents, new compounds can be prepared in SPE columns. For example, oxidation reactions can be performed by mixing the silica with the appropriate oxidizing agents. Aldol condensation reactions can also be conducted in SPE columns. In another type of application, SPE has been adopted as an alternative to liquid–liquid extraction.

² This is a Strata SPE column available from Phenomenex, 411 Madrid Ave., Torrance, CA 90501-1430; phone (310) 212-0555. Part number: 8B-S001-JCH-S, Strata C-18-E, 1000 mg sorbent/6-mL tube.

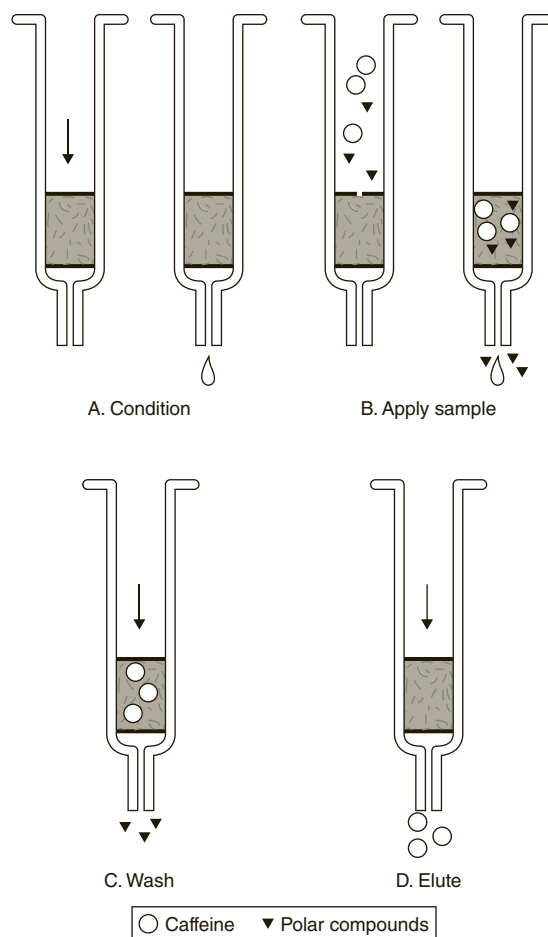



Figure 12.18
Steps to remove caffeine from tea or coffee.

PROBLEMS


- Suppose solute A has a distribution coefficient of 1.0 between water and diethyl ether. Demonstrate that if 4.0 mL of a solution of 0.20 g of A in water were extracted with two 1.0-mL portions of ether, a smaller amount of A would remain in the water than if the solution were extracted with one 2.0-mL portion of ether.
- Write an equation to show how you could recover the parent compounds from their respective salts (1, 2, and 4) shown in Figure 12.10.
- Aqueous hydrochloric acid was used *after* the sodium bicarbonate and sodium hydroxide extractions in the separation scheme shown in Figure 12.10. Is it possible to use this reagent earlier in the separation scheme to achieve the same overall result? If so, explain where you would perform this extraction.
- Using aqueous hydrochloric acid, sodium bicarbonate, or sodium hydroxide solutions, devise a separation scheme using the style shown in Figure 12.10 to separate the following two-component mixtures. All the substances are soluble in ether. Also indicate how you would recover each of the

compounds from its respective salts. Outline your procedure using a flow-chart (see Section 12.12).


- a. Give two different methods for separating this mixture.

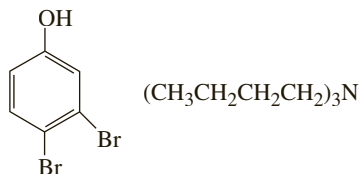

- b. Give two different methods for separating this mixture.

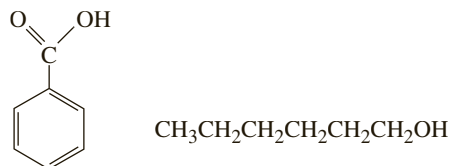

- c. Give one method for separating this mixture.

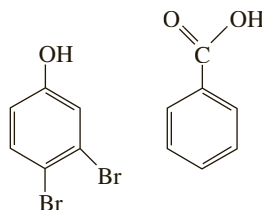

5. Solvents other than those in Table 12.1 may be used for extractions. Determine the relative positions of the organic layer and the aqueous layer in a conical vial or separatory funnel after shaking each of the following solvents with an aqueous phase. Find the densities for each of these solvents in a handbook (see Technique 4).
- 1,1,1-Trichloroethane
 - Hexane
6. A student prepares ethyl benzoate by the reaction of benzoic acid with ethanol using a sulfuric acid catalyst. The following compounds are found in the crude reaction mixture: ethyl benzoate (major component), benzoic acid, ethanol, and sulfuric acid. Using a handbook, obtain the solubility properties in water for each of these compounds (see Technique 4). Indicate how you would remove benzoic acid, ethanol, and sulfuric acid from ethyl benzoate. At some point in the purification, you should also use an aqueous sodium bicarbonate solution. Outline your procedure using a flowchart (see Section 12.12).
7. Calculate the weight of water that could be removed from a wet organic phase using 50.0 mg of magnesium sulfate. Assume that it gives the hydrate listed in Table 12.2.
8. Explain exactly what you would do when performing the following laboratory instructions; that is, write a short procedure.
- “Wash the organic layer with 1.0 mL of 1 M aqueous sodium bicarbonate.”
 - “Extract the aqueous layer three times with 1-mL portions of methylene chloride.”
9. Just prior to drying an organic layer with a drying agent, you notice water droplets in the organic layer. What should you do next?
10. What should you do if there is some question about which layer is the organic one during an extraction procedure?

11. Saturated aqueous sodium chloride ($d = 1.2 \text{ g/mL}$) is added to the following mixtures in order to dry the organic layer. Which layer is likely to be on the bottom in each case?
- Sodium chloride layer or a layer containing a high-density organic compound dissolved in methylene chloride ($d = 1.4 \text{ g/mL}$)
 - Sodium chloride layer or a layer containing a low-density organic compound dissolved in methylene chloride ($d = 1.1 \text{ g/mL}$)

Physical Constants of Liquids: The Boiling Point and Density

PART A. BOILING POINTS AND THERMOMETER CORRECTION

13.1 The Boiling Point As a liquid is heated, the vapor pressure of the liquid increases to the point at which it just equals the applied pressure (usually atmospheric pressure). At this point, the liquid is observed to boil. The normal boiling point is measured at 760 mm Hg (760 torr), or 1 atm. At a lower applied pressure, the vapor pressure needed for boiling is also lowered, and the liquid boils at a lower temperature. The relation between applied pressure and temperature of boiling for a liquid is determined by its vapor pressure–temperature behavior. Figure 13.1 is an idealization of the typical vapor pressure–temperature behavior of a liquid.

Because the boiling point is sensitive to pressure, it is important to record the barometric pressure when determining a boiling point if the determination is being conducted at an elevation significantly above or below sea level. Normal atmospheric variations may affect the boiling point, but they are usually of minor importance. However, if a boiling point is being monitored during the course of a vacuum distillation (Technique 16) that is being performed with an aspirator or a vacuum pump, the variation from the atmospheric value will be especially marked. In these cases, it is important to know the pressure as accurately as possible.

As a rule of thumb, the boiling point of many liquids drops about 0.5°C for a 10-mm decrease in pressure when in the vicinity of 760 mm Hg. At lower pressures, a 10°C drop in boiling point is observed for each halving of the pressure. For example, if the observed boiling point of a liquid is 150°C at 10 mm pressure, then the boiling point would be about 140°C at 5 mm Hg.

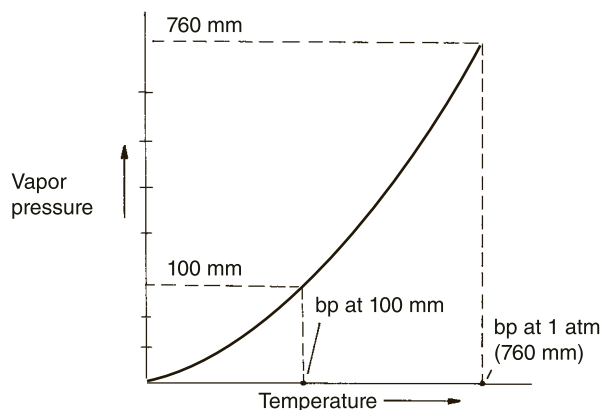


Figure 13.1
The vapor pressure–temperature curve for a typical liquid.

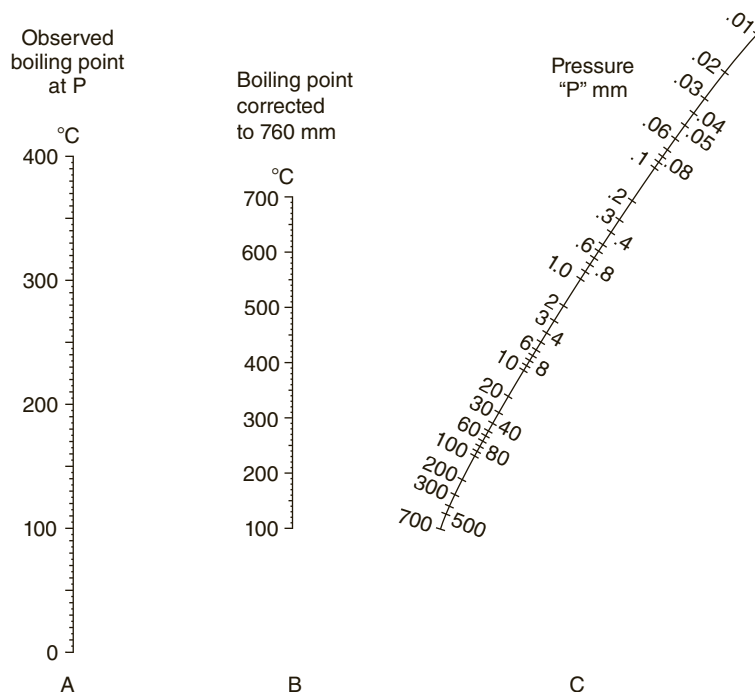


Figure 13.2

Pressure–temperature alignment nomograph. How to use the nomograph: Assume a reported boiling point of 100°C (column A) at 1 mm. To determine the boiling point at 18 mm, connect 100°C (column A) to 1 mm (column C) with a transparent plastic rule and observe where this line intersects column B (about 280°C). This value would correspond to the normal boiling point. Next, connect 280°C (column B) with 18 mm (column C) and observe where this intersects column A (151°C). The approximate boiling point will be 151°C at 18 mm. (Reprinted courtesy of EMD Chemicals, Inc.)

A more accurate estimate of the change in boiling point with a change of pressure can be made by using a nomograph. In Figure 13.2, a nomograph is given, and a method is described for using it to obtain boiling points at various pressures when the boiling point is known at some other pressure.

13.2 Determining the Boiling Point—Microscale and Macroscale Methods

Several experimental methods of determining boiling points are available. Selecting a method depends on how much liquid is available and the availability of specific apparatus. In either microscale or macroscale experiments in which 0.3–0.5 mL of liquid is available, the semimicroscale direct method is usually most reliable. If less material is available, it will be necessary to perform either the semimicroscale or microscale inverted capillary method. With practice, these methods can be reliable, too. With larger quantities in both microscale and macroscale experiments, the boiling point can also be observed while performing a distillation.

Semimicroscale Direct Method. The apparatus for this method is shown in Figure 13.3. With this method, the bulb of the thermometer can be immersed in vapor from the boiling liquid for a period long enough to allow it to equilibrate and give a good temperature reading. A 13-mm × 100-mm test tube works well in this procedure. Use 0.5 mL of liquid and a small, inert carborundum (black) boiling

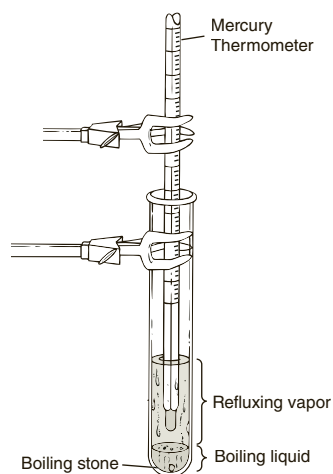


Figure 13.3
Macroscale method of determining the boiling point.

stone. This method works best with a partial immersion (76 mm) mercury thermometer (see Section 13.3). It is not necessary to perform a stem correction with this type of thermometer. This method also works well with a digital thermometer (see Section 13.4).

Place the bulb of the thermometer as close as possible to the boiling liquid without actually touching it. The best heating device is a hot plate with either an aluminum block or a sand bath.¹

While you are heating the liquid, it is helpful to record the temperature at 1-minute intervals. This makes it easier to keep track of changes in the temperature and to know when the liquid has reached the boiling point. The liquid must boil vigorously, such that you see a reflux ring above the bulb of the thermometer and drops of liquid condensing on the sides of the test tube. Note that with some liquids, the reflux ring will be faint, and you must look closely to see it. The boiling point is reached when the temperature reading on the thermometer has remained constant at its highest observed value for 2–3 minutes. It is usually best to turn the heat control on the hot plate to a relatively high setting initially, especially if you are starting with a cold hot plate and aluminum block or sand bath. If the temperature begins to level off at a relatively low temperature (less than about

100°C) or if the reflux ring reaches the immersion ring on the thermometer, you should turn down the heat-control setting immediately.

Two problems can occur when you perform this boiling-point procedure. The first is much more common and occurs when the temperature appears to be leveling off at a temperature below the boiling point of the liquid. This is more likely to happen with a relatively high-boiling liquid (boiling points greater than about 150°C) or when the sample is not heated sufficiently. The best way to prevent this problem is to heat the sample more strongly. With high-boiling liquids, it may be helpful to wait for the temperature to remain constant for 3–4 minutes to make sure that you have reached the actual boiling point.

The second problem, which is rare, occurs when the liquid evaporates completely, and the temperature inside the dry test tube may rise higher than the actual boiling point of the liquid. This is more likely to happen with low-boiling liquids (boiling point less than 100°C) or if the temperature on the hot plate is set too high for too long. To check for this possibility, observe the amount of liquid remaining in the test tube as soon as you have finished with the procedure. If there is no liquid remaining, it is possible that the highest temperature you observed is greater than the boiling point of the liquid. In this case, you should repeat the boiling-point determination, heating the sample less strongly or using more sample.

Depending on the skill of the person performing this technique, boiling points may sometimes be inaccurate. When experimental boiling points are inaccurate, it is more common for them to be lower than the literature value, and inaccuracies are more likely to occur for higher-boiling liquids. Carefully following the previous instructions will make it more likely that your experimental value will be close to the literature value.

With smaller amounts of material, you can carry out a microscale or semimicroscale determination of the boiling point by using the apparatus shown in Figure 13.4.

¹ *Note to the Instructor:* The aluminum block should have a hole drilled in it that goes all the way through the block and is just slightly larger than the outside diameter of the test tube. A sand bath can be conveniently prepared by adding 40 mL of sand to a 150-mL beaker or by using a heating mantle partially filled with sand. For additional comments about these heating methods, see the Instructor's Manual, Experiment 7, "Infrared Spectroscopy and Boiling-Point Determination."

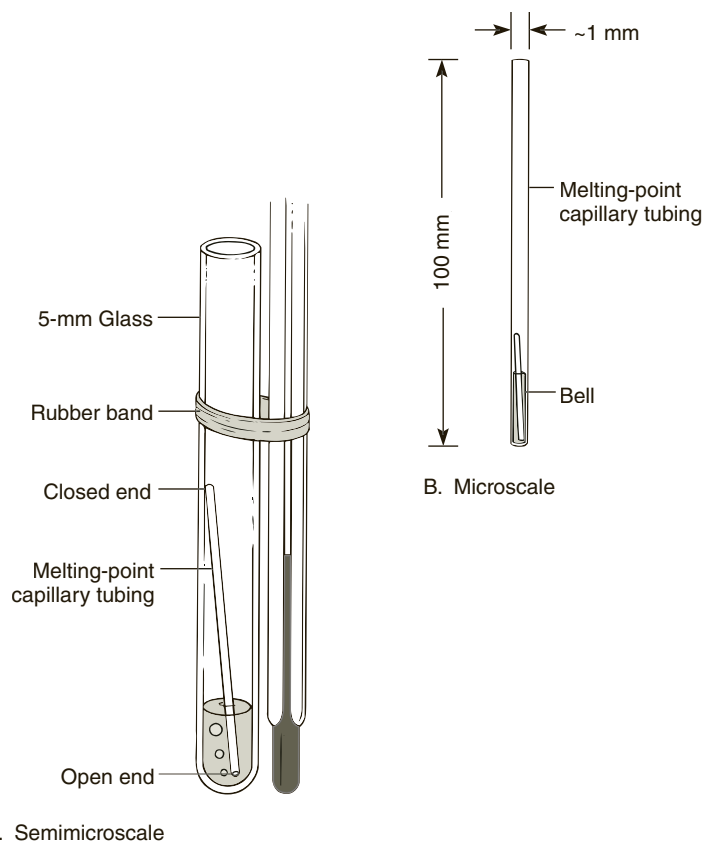


Figure 13.4
Boiling-point determinations.

Semimicroscale Inverted Capillary Method. To carry out the semimicroscale determination, attach a piece of 5-mm glass tubing (sealed at one end) to a thermometer with a rubber band or a thin slice of rubber tubing. The liquid whose boiling point is being determined is introduced with a Pasteur pipette into this piece of tubing, and a short piece of melting-point capillary (sealed at one end) is dropped in with the open end down. The whole unit is then placed in a Thiele tube. The rubber band should be placed above the level of the oil in the Thiele tube; otherwise, the band may soften in the hot oil. When positioning the band, keep in mind that the oil will expand when heated. Next, the Thiele tube is heated in the same fashion as described in Technique 9, Section 9.6, for determining a melting point. Heating is continued until a rapid and continuous stream of bubbles emerges from the inverted capillary. At this point, you should stop heating. Soon, the stream of bubbles slows down and stops. When the bubbles stop, the liquid enters the capillary tube. The moment at which the liquid enters the capillary tube corresponds to the boiling point of the liquid, and the temperature is recorded.

Microscale Inverted Capillary Method. In microscale experiments, there often is too little product available to use the semimicroscale method just described. However, the method can be scaled down in the following manner. The liquid is placed in a 1-mm melting-point capillary tube to a depth of about 4–6 mm (see Figure 13.4B). Use a syringe or a Pasteur pipette that has had its tip drawn thinner to transfer the liquid into the capillary tube. It may be necessary to use a centrifuge to transfer the liquid to the bottom of the tube. Next, prepare an appropriately sized inverted capillary, or **bell**.

The easiest way to prepare a bell is to use a commercial micropipette, such as a 10- μL Drummond “microcap.” These are available in vials of 50 or 100 microcaps and are inexpensive. To prepare the bell, cut the microcap in half with a file or scorer and then seal one end by inserting it a small distance into a flame, turning it on its axis until the opening closes.

If microcaps are not available, a piece of 1-mm open-end capillary tubing (same size as a melting-point capillary) can be rotated along its axis in a flame while being held horizontally. Use your index fingers and thumbs to rotate the tube; do not change the distance between your two hands while rotating. When the tubing is soft, remove it from the flame and pull it to a thinner diameter. When pulling, keep the tube straight by *moving both your hands and your elbows outward* by about 4 inches. Hold the pulled tube in place a few moments until it cools. Using the edge of a file or your fingernail, break out the thin center section. Seal one end of the thin section in the flame; then break it to a length that is about one and one-half times the height of your sample liquid (6–9 mm). Be sure the break is done squarely. Invert the bell (open end down), and place it in the capillary tube containing the sample liquid. Push the bell to the bottom with a fine copper wire if it adheres to the side of the capillary tube. A centrifuge may be used if you prefer. Figure 13.5 shows the construction method for the bell and the final assembly.

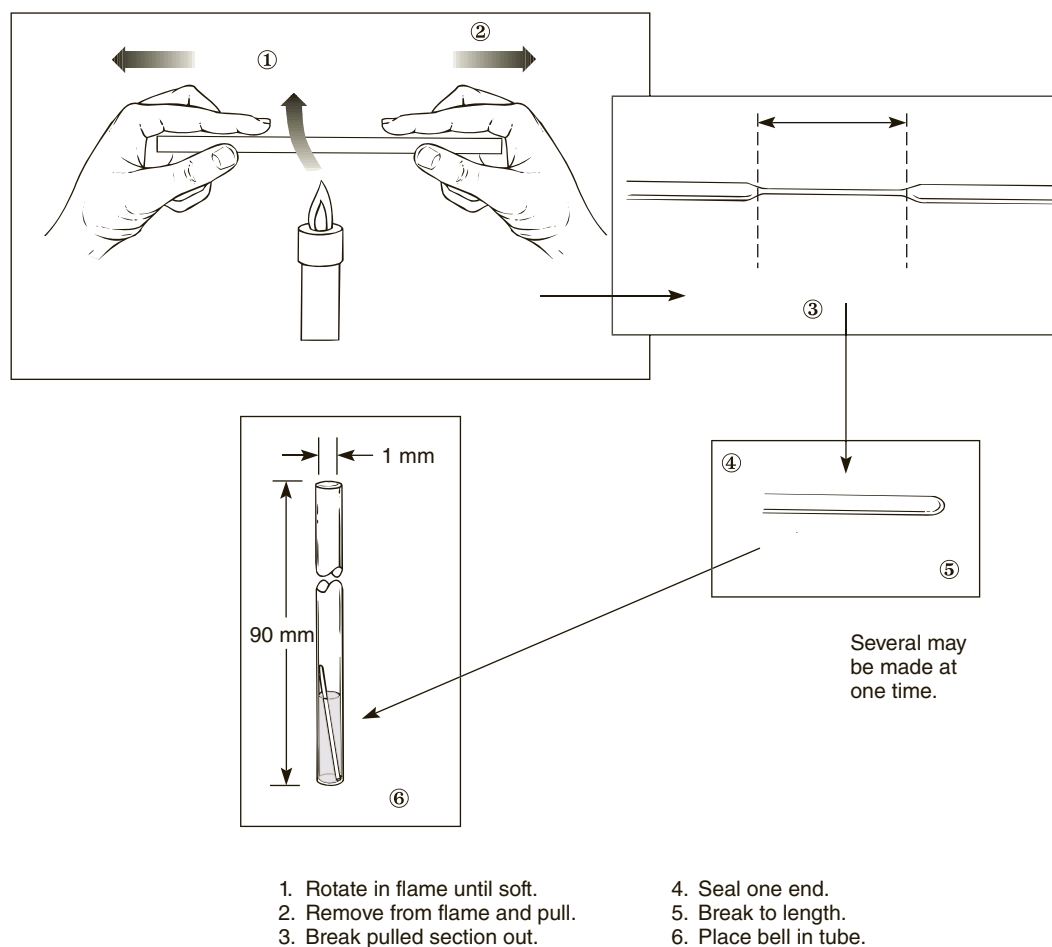


Figure 13.5
Construction of a microcapillary bell for microscale boiling-point determination.

Place the microscale assembly in a standard melting-point apparatus (or a Thiele tube if an electrical apparatus is not available) to determine the boiling point. Heating is continued until a rapid and continuous stream of bubbles emerges from the inverted capillary. At this point, stop heating. Soon, the stream of bubbles slows down and stops. When the bubbles stop, the liquid enters the capillary tube. The moment at which the liquid enters the capillary tube corresponds to the boiling point of the liquid, and the temperature is recorded.

Explanation of the Method. During the initial heating, the air trapped in the inverted bell expands and leaves the tube, giving rise to a stream of bubbles. When the liquid begins boiling, most of the air has been expelled; the bubbles of gas are due to the boiling action of the liquid. Once the heating is stopped, most of the vapor pressure left in the bell comes from the vapor of the heated liquid that seals its open end. There is always vapor in equilibrium with a heated liquid. If the temperature of the liquid is above its boiling point, the pressure of the trapped vapor will either exceed or equal the atmospheric pressure. As the liquid cools, its vapor pressure decreases. When the vapor pressure drops just below atmospheric pressure (just below the boiling point), the liquid is forced into the capillary tube.

Difficulties. Three problems are common to this method. The first arises when the liquid is heated so strongly that it evaporates or boils away. The second arises when the liquid is not heated above its boiling point before heating is discontinued. If the heating is stopped at any point below the actual boiling point of the sample, the liquid enters the bell *immediately*, giving an apparent boiling point that is too low. Be sure you observe a continuous stream of bubbles, too fast for individual bubbles to be distinguished, before lowering the temperature. Also be sure the bubbling action decreases slowly before the liquid enters the bell. If your melting-point apparatus has fine enough control and fast response, you can actually begin heating again and force the liquid out of the bell before it becomes completely filled with the liquid. This allows a second determination to be performed on the same sample. The third problem is that the bell may be so light that the bubbling action of the liquid causes the bell to move up the capillary tube. This problem can sometimes be solved by using a longer (heavier) bell or by sealing the bell so that a larger section of solid glass is formed at the sealed end of the bell.

When measuring temperatures above 150°C, thermometer errors can become significant. For an accurate boiling point with a high-boiling liquid, you may wish to apply a **stem correction** to the thermometer, as described in Section 13.3, or to calibrate the thermometer, as described in Technique 9, Section 9.9.

Microscale or Macroscale–Distillation Method. When you have large quantities of material, you can simply record the boiling point (or boiling range) as viewed on a thermometer while performing a simple distillation (see Technique 14). When this method is used to determine a boiling point, it is best to use a partial immersion mercury thermometer or a digital thermometer for more accurate readings.

13.3 Glass Thermometers and Stem Corrections

Three types of glass thermometers are available: bulb immersion, partial immersion (stem immersion), and total immersion. Bulb immersion thermometers are calibrated by the manufacturer to give correct temperature readings when only the bulb (not the rest of the thermometer) is placed in the medium to be measured. Partial immersion thermometers are calibrated to give correct temperature readings when they are immersed to a specified depth in the medium to be measured. Partial immersion thermometers are easily recognized because the manufacturer always

scores a mark, or immersion ring, completely around the stem at the specified depth of immersion. The immersion ring is normally found below any of the temperature calibrations. Total immersion thermometers are calibrated when the entire thermometer is immersed in the medium to be measured. The three types of thermometer are often marked on the back (opposite side from the calibrations) by the words *bulb*, *immersion*, or *total*, but this may vary from one manufacturer to another.

Boiling-point determination and distillation are two techniques in which an accurate temperature reading may be obtained most easily with a partial immersion thermometer. A common immersion length for this type of thermometer is 76 mm. This length works well for these two techniques because the hot vapors are likely to surround the bottom of the thermometer up to a point fairly close to the immersion line. If a total immersion thermometer is used in these applications, a stem correction, which is described later, must be used to obtain an accurate temperature reading.

The liquid used in thermometers may be either mercury or a colored organic liquid such as an alcohol. Because mercury is highly poisonous and is difficult to clean up completely when a thermometer is broken, many laboratories now use non-mercury thermometers. When a highly accurate temperature reading is required, such as in a boiling-point determination or in some distillations, mercury thermometers may have an advantage over non-mercury thermometers for two reasons. Mercury has a lower coefficient of expansion than the liquids used in non-mercury thermometers. Therefore, a partial immersion mercury thermometer will give a more accurate reading when the thermometer is not immersed in the hot vapors exactly to the immersion line. In other words, the mercury thermometer is more forgiving. Furthermore, because mercury is a better conductor of heat, a mercury thermometer will respond more quickly to changes in the temperature of the hot vapors. If the temperature is read before the thermometer reading has stabilized, which is more likely to occur with a non-mercury thermometer, the temperature reading will be inaccurate.

Manufacturers design total immersion thermometers to read correctly only when they are immersed totally in the medium to be measured. The entire mercury thread must be covered. Because this situation is rare, a **stem correction** should be added to the observed temperature. This correction, which is positive, can be fairly large when high temperatures are being measured. Keep in mind, however, that if your thermometer has been calibrated for its desired use (such as described in Technique 9, Section 9.9 for a melting-point apparatus), a stem correction should not be necessary for any temperature within the calibration limits. You are most likely to want a stem correction when you are performing a distillation. If you determine a melting point or boiling point using an uncalibrated, total immersion thermometer, you will also want to use a stem correction.

When you wish to make a stem correction for a total immersion thermometer, the following formula may be used. It is based on the fact that the portion of the mercury thread in the stem is cooler than the portion immersed in the vapor or the heated area around the thermometer. The mercury will not have expanded in the cool stem to the same extent as in the warmed section of the thermometer. The equation used is

$$(0.000154)(T - t_1)(T - t_2) = \text{correction to be added to } T \text{ observed}$$

1. The factor 0.000154 is a constant, the coefficient of expansion for the mercury in the thermometer.
2. The term $T - t_1$ corresponds to the length of the mercury thread not immersed in the heated area. Use the temperature scale on the thermometer itself for this measurement rather than an actual length unit. T is the observed temperature, and t_1 is the *approximate* place where the heated part of the stem ends and the cooler part begins.

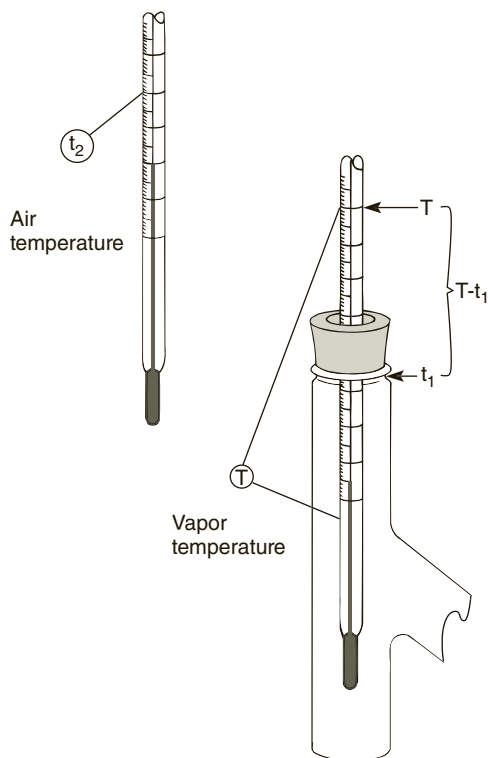


Figure 13.6
Measurement of a thermometer stem correction during distillation.

3. The term $T - t_2$ corresponds to the difference between the temperature of the mercury in the vapor T and the temperature of the mercury in the air outside the heated area (room temperature). The term T is the observed temperature, and t_2 is measured by hanging another thermometer so the bulb is close to the stem of the main thermometer.

Figure 13.6 shows how to apply this method for a distillation. By the formula just given, it can be shown that high temperatures are more likely to require a stem correction and that low temperatures need not be corrected. The following sample calculations illustrate this point.

Example 1	Example 2
$T = 200^\circ\text{C}$	$T = 100^\circ\text{C}$
$t_1 = 0^\circ\text{C}$	$t_1 = 0^\circ\text{C}$
$t_2 = 35^\circ\text{C}$	$t_2 = 35^\circ\text{C}$
$(0.000154)(200)(165) = 5.1^\circ$ stem correction	$(0.000154)(100)(65) = 1.0^\circ$ stem correction
$200^\circ\text{C} + 5^\circ\text{C} = 205^\circ\text{C}$ corrected temperature	$100^\circ\text{C} + 1^\circ\text{C} = 101^\circ\text{C}$ corrected temperature

13.4 Digital Thermometers

Rather than using a glass thermometer to determine a boiling point or to monitor the temperature during a distillation, one can use a digital thermometer with a stainless steel temperature probe (see Figure 13.7). When a digital thermometer is used with the semimicroscale direct method of determining a boiling point (see Semimicroscale Direct Method in Section 13.2), the boiling point can usually be

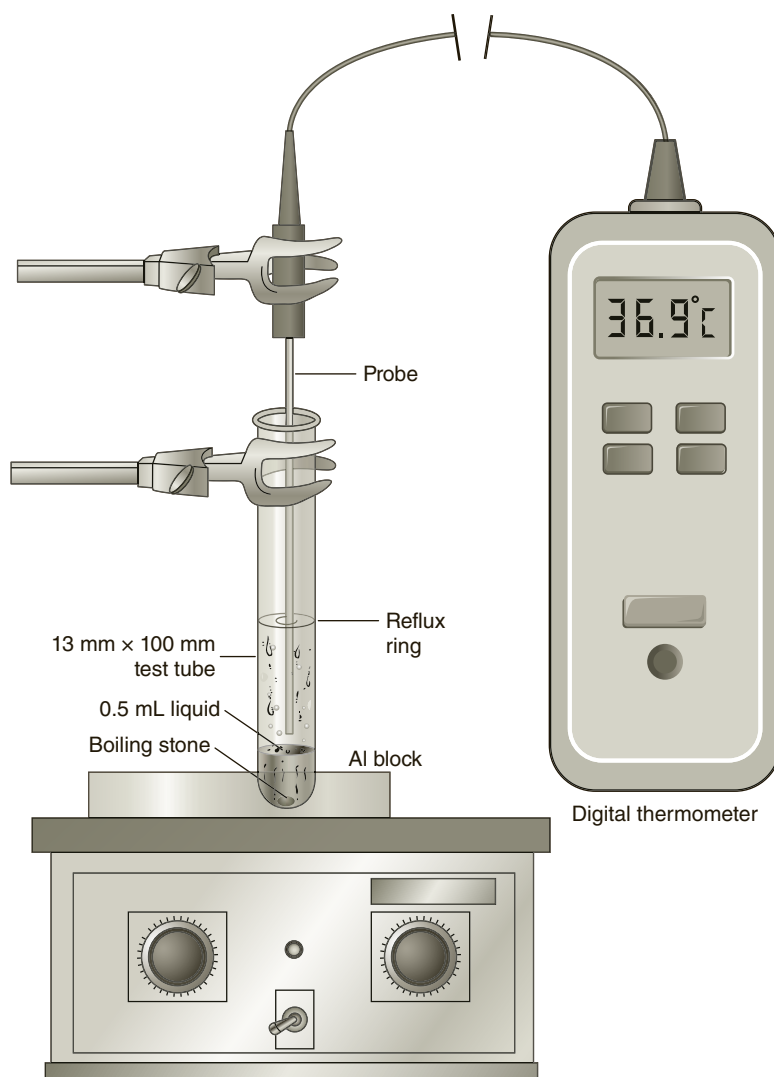


Figure 13.7
Boiling point determination with a digital thermometer.

determined to within 2-3 degrees of the literature value. The temperature probe (or thermocouple) works only in a given temperature range. It is therefore important to select a probe that has a maximum temperature that is higher than the boiling points of the liquids that you will be measuring. See the Instructor's Manual, Experiment 8, Simple and Fractional Distillation, for more specific information about selecting an appropriate temperature probe.

An even better way to measure the temperature is to use the Vernier LabQuest unit, with a stainless steel temperature probe. The LabQuest is a handheld device that can be used as a stand-alone unit. As with the digital thermometers, the temperature probe used with these devices must also work in a range that exceeds the boiling points of the liquids you will be attempting to measure. The data (temperature vs. time) are displayed either on the screen of the LabQuest unit while it is being collected.

Being able to see a graph of temperature vs. time when performing a distillation gives you a better sense of when the different liquids are distilling and makes it easier to control the distillation (see Technique 14, Section 14.5).

PART B. DENSITY

13.5 Density

Density is defined as mass per unit volume and is generally expressed in units of grams per milliliter (g/mL) for a liquid and grams per cubic centimeter (g/cm³) for a solid.

$$\text{Density} = \frac{\text{mass}}{\text{volume}} \quad \text{or} \quad D = \frac{M}{V}$$

In organic chemistry, density is most commonly used in converting the weight of liquid to a corresponding volume, or vice versa. It is often easier to measure a volume of a liquid than to weigh it. As a physical property, density is also useful for identifying liquids in much the same way that boiling points are used.

Although precise methods that allow the measurements of the densities of liquids at the microscale level have been developed, they are often difficult to perform. An approximate method for measuring densities can be found in using a 100- μL (0.100-mL) automatic pipette (Technique 5, Section 5.1). Clean, dry, and preweigh one or more conical vials (including their caps and liners) and record their weights. Handle these vials with a tissue to avoid getting your fingerprints on them. Adjust the automatic pipette to deliver 100 μL and fit it with a clean, new tip. Use the pipette to deliver 100 μL of the unknown liquid to each of your tared vials. Cap them so that the liquid does not evaporate. Reweigh the vials and use the weight of the 100 μL of liquid delivered to calculate a density for each case. It is recommended that from three to five determinations be performed, that the calculations be performed to three significant figures, and that all the calculations be averaged to obtain the final result. This determination of the density will be accurate to within two significant figures. Table 13.1 compares some literature values with those that could be obtained by this method.

Table 13.1 Densities determined by the automatic pipette method (g/mL)

<i>Substance</i>	<i>BP</i>	<i>Literature</i>	<i>100 μL</i>
Water	100	1.000	1.01
Hexane	69	0.660	0.66
Acetone	56	0.788	0.77
Dichloromethane	40	1.330	1.27
Diethyl ether	35	0.713	0.67

PROBLEMS

- Using the pressure–temperature alignment chart in Figure 13.2, answer the following questions.
 - What is the normal boiling point (at 760 mm Hg) for a compound that boils at 150°C at 10 mm Hg pressure?
 - At what temperature would the compound in (a) boil if the pressure were 40 mm Hg?
 - A compound was distilled at atmospheric pressure and had a boiling point of 285°C. What would be the approximate boiling range for this compound at 15 mm Hg?

2. Calculate the corrected boiling point for nitrobenzene by using the method given in Section 13.3. The boiling point was determined using an apparatus similar to that shown in Figure 13.3. Assume that a total immersion thermometer was used. The observed boiling point was 205°C. The reflux ring in the test tube just reached up to the 0°C mark on the thermometer. A second thermometer suspended alongside the test tube, at a slightly higher level than the one inside, gave a reading of 35°C.
3. Suppose that you had calibrated the thermometer in your melting-point apparatus against a series of melting-point standards. After reading the temperature and converting it using the calibration chart, should you also apply a stem correction? Explain.
4. The density of a liquid was determined by the automatic pipette method. A 100- μ L automatic pipette was used. The liquid had a mass of 0.082 g. What was the density in grams per milliliter of the liquid?
5. During the microscale boiling-point determination (see Section 13.2, Microscale Inverted Capillary Method) of an unknown liquid, heating was discontinued at 154°C and the liquid immediately began to enter the inverted bell. Heating was begun again at once, and the liquid was forced out of the bell. Heating was again discontinued at 165°C, at which time a rapid stream of bubbles emerged from the bell. On cooling, the rate of bubbling gradually diminished until the liquid reached a temperature of 161°C and entered and filled the bell. Explain this sequence of events. What was the boiling point of the liquid?

Simple Distillation

Distillation is the process of vaporizing a liquid, condensing the vapor, and collecting the condensate in another container. This technique is useful for separating a liquid mixture when the components have different boiling points or when one of the components will not distill. It is one of the principal methods of purifying a liquid. Four basic distillation methods are available to the chemist: simple distillation, vacuum distillation (distillation at reduced pressure), fractional distillation, and steam distillation. This technique chapter will discuss simple distillation. Vacuum distillation will be discussed in Technique 16. Fractional distillation will be discussed in Technique 15, and steam distillation will be discussed in Technique 18.

14.1 The Evolution of Distillation Equipment

There are probably more types and styles of distillation apparatus than exist for any other technique in chemistry. Over the centuries, chemists have devised just about every conceivable design. The earliest known types of distillation apparatus were the **alembic** and the **retort** (Figure 14.1). They were used by alchemists in the Middle Ages and the Renaissance and probably even earlier by Arabic chemists. Most other distillation equipment has evolved as variations on these designs.

Figure 14.1 shows several stages in the evolution of distillation equipment as it relates to the organic laboratory. It is not intended to be a complete history; rather, it is representative. Up until recent years, equipment based on the retort design was common in the laboratory. Although the retort itself was still in use early in the twentieth century, it had evolved by that time into the distillation flask and water-cooled condenser combination. This early equipment was connected with corks. By 1958, most introductory laboratories were beginning to use “organic lab kits” that included glassware connected by standard-taper glass joints. The original lab kits contained large $\text{F } 24/40$ joints. Within a short time, they became smaller, with $\text{F } 19/22$ and even $\text{F } 14/20$ joints. These later kits are still being used today in many organic courses. Small-scale variations of these kits are also used today by chemical researchers, but they are too expensive to use in an introductory laboratory. Instead, the “microscale” equipment you are using in this course is coming into common use. This equipment has $\text{F } 14/10$ standard-taper joints, threaded outer joints with screwcap connectors, and an internal O-ring. The distillation apparatus in microscale kits is designed for work with small amounts of material, and it is different from the more traditional larger-scale equipment. It is perhaps more closely related to the alembic design than to that of the retort. Because both types of equipment are in use today, after we describe microscale equipment, we will also show the equivalent large-scale apparatus used to perform distillation.

14.2 Distillation Theory

In the traditional distillation of a pure substance, vapor rises from the distillation flask and comes into contact with a thermometer that records its temperature. The vapor then passes through a condenser, which relieves the vapor and passes it into the receiving flask. The temperature observed during the distillation of a **pure**

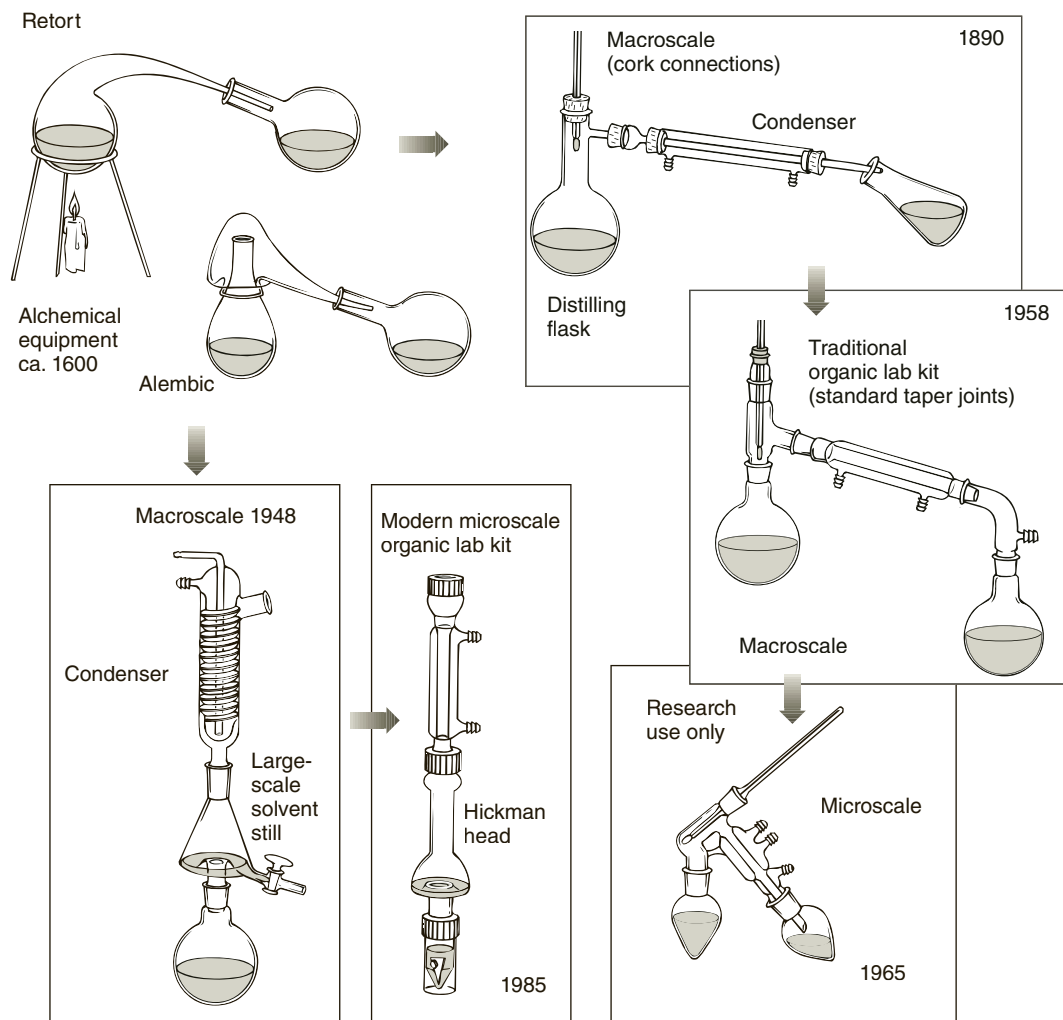


Figure 14.1 Some stages in the evolution of microscale distillation equipment from alchemical equipment (dates represent approximate time of popular use).

substance remains constant throughout the distillation so long as both vapor *and* liquid are present in the system (see Figure 14.2A). When a **liquid mixture** is distilled, often the temperature does not remain constant but increases throughout the distillation. The reason for this is that the composition of the vapor that is distilling varies continuously during the distillation (see Figure 14.2B).

For a liquid mixture, the composition of the vapor in equilibrium with the heated solution is different from the composition of the solution itself. This is shown in Figure 14.3, which is a phase diagram of the typical vapor–liquid relationship for a two-component system (A + B).

On this diagram, horizontal lines represent constant temperatures. The upper curve represents vapor composition, and the lower curve represents liquid composition. For any horizontal line (constant temperature), like that shown at t , the intersections of the line with the curves give the compositions of the liquid and the vapor that are in equilibrium with each other at that temperature. In the diagram, at temperature t , the intersection of the curve at X indicates that liquid of composition W will be in equilibrium with vapor of composition Z, which corresponds to the intersection at Y. Composition is given as a mole percentage of A and B in the mixture. Pure A, which boils at temperature t_A , is represented at the left. Pure B, which boils

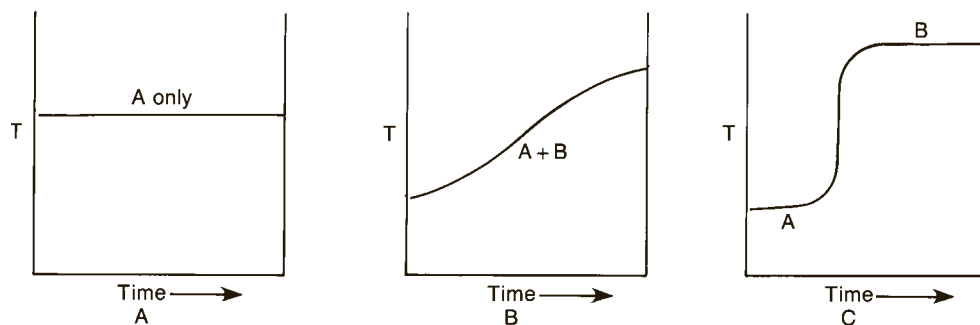


Figure 14.2

Three types of temperature behavior during a simple distillation. A. A single pure component. B. Two components of similar boiling points. C. Two components with widely differing boiling points. Good separations are achieved in A and C.

at temperature t_B , is represented on the right. For either pure A or pure B, the vapor and liquid curves meet at the boiling point. Thus, either pure A or pure B will distill at a constant temperature (t_A or t_B). Both the vapor and the liquid must have the same composition in either of these cases. This is not the case for mixtures of A and B.

A mixture of A and B of composition W will have the following behavior when heated. The temperature of the liquid mixture will increase until the boiling point of the mixture is reached. This corresponds to following line WX from W to X, the boiling point of the mixture t . At temperature t the liquid begins to vaporize, which corresponds to line XY. The vapor has the composition corresponding to Z. In other words, the first vapor obtained in distilling a mixture of A and B does not consist of pure A. It is richer in A than the original mixture but still contains a significant amount of the higher-boiling component B, even from the very beginning of the distillation. The result is that it is never possible to separate a mixture completely by a simple distillation. However, in two cases it is possible to get an acceptable separation into relatively pure components. In the first case, if the boiling points of A and B differ by a large amount (>100 degrees) and if the distillation is carried out carefully, it will be possible to get a fair separation of A and B. In the second case, if A contains a fairly small amount of B ($<10\%$), a reasonable separation of A from B can be achieved. When the boiling-point

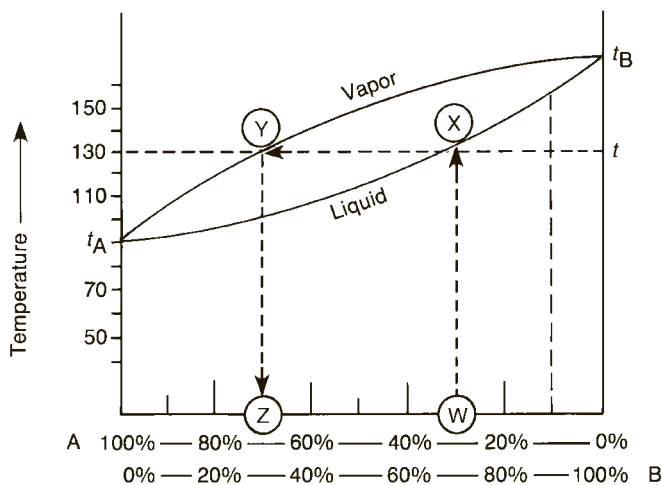


Figure 14.3
Phase diagram for a typical liquid mixture of two components.

differences are not large and when highly pure components are desired, it is necessary to perform a **fractional distillation**. Fractional distillation is described in Technique 15, in which the behavior during a simple distillation is also considered in detail. Note only that as vapor distills from the mixture of composition W (Figure 14.3), it is richer in A than is the solution. Thus, the composition of the material left behind in the distillation becomes richer in B (moves to the right from W toward pure B in the graph). A mixture of 90% B (dotted line on the right side in Figure 14.3) has a higher boiling point than at W. Hence, the temperature of the liquid in the distillation flask will increase during the distillation, and the composition of the distillate will change (as is shown in Figure 14.3B).

When two components that have a large boiling-point difference are distilled, the temperature

remains constant while the first component distills. If the temperature remains constant, a relatively pure substance is being distilled. After the first substance distills, the temperature of the vapors rises, and the second component distills, again at a constant temperature. This is shown in Figure 14.2C. A typical application of this type of distillation might be an instance of a reaction mixture containing the desired component A (bp 140°C) contaminated with a small amount of undesired component B (bp 250°C) and mixed with a solvent such as diethyl ether (bp 36°C). The ether is removed easily at low temperature. Pure A is removed at a higher temperature and collected in a separate receiver. Component B can then be distilled, but it usually is left as a residue and not distilled. This separation is not difficult and represents a case where simple distillation might be used to advantage.

14.3 Microscale Equipment

Most large-scale distillation equipment requires the distilled liquid to travel a long distance from the distillation flask, through the condenser, to the receiving flask. When working at the microscale level, a long distillation path must be avoided. With small quantities of liquid, there are too many opportunities to lose all the sample. The liquid will adhere to, or wet, surfaces and get lost in every little nook and cranny of the system. A system with a long path also has a large volume, and a small amount of liquid may not produce enough vapor to fill it. Small-scale distillation requires a “short path” distillation. In order to make the distilling path as short as possible, the **Hickman head** has been adopted as the principal receiving device for most microscale distilling operations.

The Hickman Head. Two types of Hickman head (also called a Hickman “still”) are shown in Figure 14.4. One of these variations has a convenient opening, or port, in the side, making removal of liquid that has collected in it easier. In operation, the liquid to be distilled is placed in a flask or vial attached to the bottom joint of the Hickman head and heated. If desired, you can attach a condenser to the top joint. Either a magnetic spin vane or a boiling stone is used to prevent bumping. Some typical assemblies are shown in Figures 14.5 and 14.7. The vapors of the heated liquid rise upward and are cooled and condensed on either the walls of the condenser or, if no condenser is used, on the inside walls of the Hickman head itself. As liquid drains downward, it collects in the circular well at the bottom of the still.

Collecting Fractions. The liquid that distills is called the **distillate**. Portions of the distillate collected during the course of a distillation are called **fractions**. A small fraction (usually discarded) collected before the distillation is begun in earnest is called a **forerun**. The well in a Hickman head can contain anywhere from 1 to 2 mL of liquid. In the style with the side port, fractions may be removed by opening the port and inserting a Pasteur pipette (Figure 14.6C). The unported head works equally well, but the head is emptied from the top by using a Pasteur pipette (Figure 14.6A). If a condenser or an internal thermometer is used, the distilling apparatus must be partially disassembled to remove liquid when the well fills. In some stills, the inner diameter of the head is small, and it is difficult to reach in at an angle with the pipette and make contact with the liquid. To remedy this problem, you may be able to use the longer (9-inch) Pasteur pipette instead of the shorter (5¾-inch) one. The longer pipette has a much longer narrow section (tip) and can adapt more effectively to the required angle. The disadvantage of the longer tip is that you are more likely to break it off inside the still. You may prefer to modify a short pipette by bending its tip slightly in a flame (Figure 14.6B).

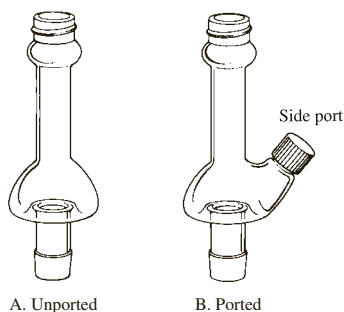


Figure 14.4
The Hickman head.

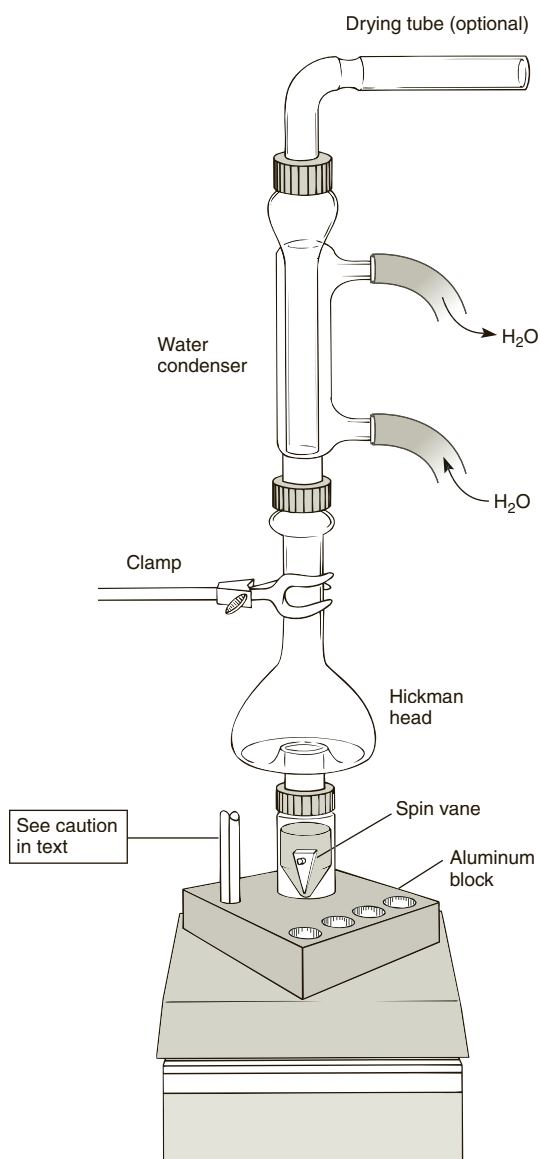


Figure 14.5
Basic microscale distillation (external monitoring of temperature). Do not use a mercury thermometer.

Choice of Condenser. If you are careful (slow heating) or if the liquid to be distilled has a high boiling point, it may not be necessary to use a condenser with the Hickman head (Figure 14.7). In this case, the liquid being distilled must condense on the cooler sides of the head itself without any being lost through evaporation. If the liquid has a low boiling point or is very volatile, a condenser must be used. With very volatile liquids, a watercooled condenser must be used; however, an air-cooled condenser may suffice for less demanding cases. When using a water condenser, remember that water should enter the lower opening and exit from the upper one. If the hoses carrying the water in and out are connected in reverse fashion, the water jacket of the condenser will not fill completely.

Sealed Systems. Whenever you perform a distillation, be sure the system you are heating is not sealed off completely from the outside atmosphere. During a distillation, the air and vapors inside the system will both expand and contract. If pressure builds up inside a sealed system, the apparatus may explode. In performing a distillation, you should leave a small opening at the far end of the system. If water vapor could be harmful to the substances being distilled, a calcium chloride drying tube may be used to protect the system from moisture. Carefully examine each system discussed to see how an opening to the outside is provided.

External Monitoring of Temperature. The simple assembly using the Hickman head shown in Figure 14.5 does not monitor the temperature inside the apparatus. Instead, the temperature is monitored externally with a thermometer placed in an aluminum block.

CAUTION



You should not use a mercury thermometer with an aluminum block. If it breaks, the mercury will vaporize on the hot surface. Instead, use a non-mercury glass thermometer, a metal dial thermometer, or a digital electronic temperature measuring device.

External monitoring of the temperature has the disadvantage that the exact temperature at which liquid distills is never known. In many cases, this does not matter or is unavoidable, and the boiling point of the distilled liquid can be checked later by performing a microboiling-point determination (Technique 13, Section 13.2).

As a rule, there is at least a 15-degree difference in temperature between the temperature of the aluminum block or sand bath and that of the liquid in the heated distillation vial or flask. However, the magnitude of this difference cannot be relied on. Keep in mind that the liquid in the vial or flask may be at a different temperature than the vapor that is distilling. In many procedures in this text, the *approximate* temperature of the heating device will be given instead of the boiling point of the liquid involved. Because this method of monitoring the temperature is rather approximate, you will need to make the actual heater setting based on what is supposed to be occurring in the vial or flask.

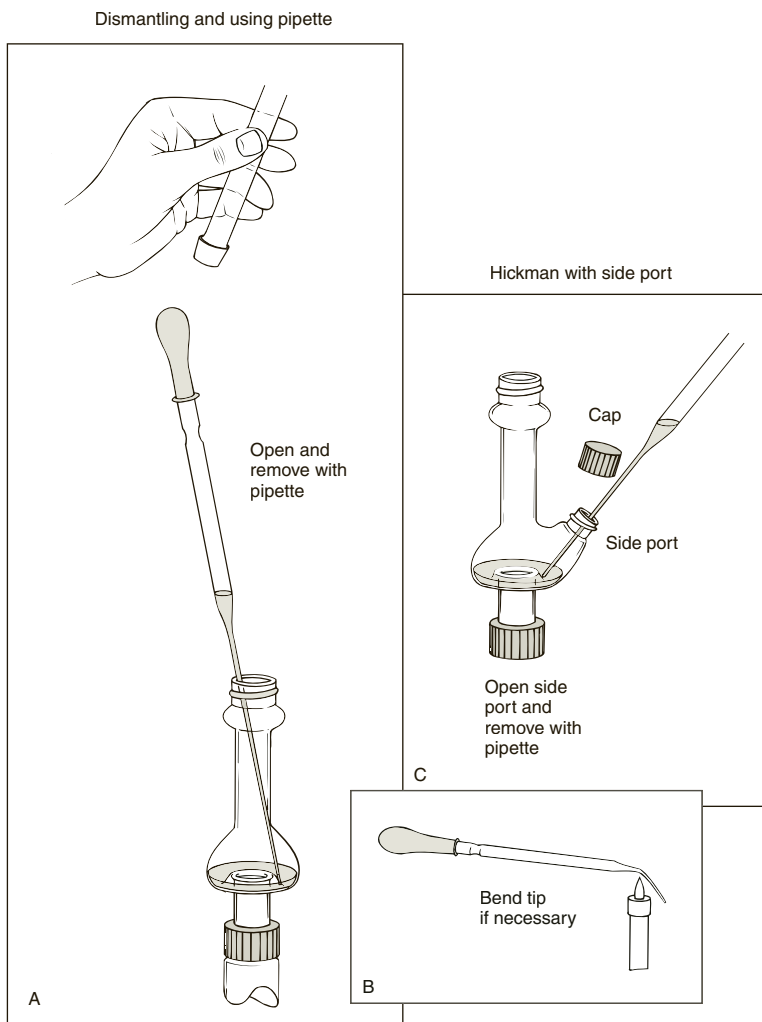


Figure 14.6
Removing fractions.

Internal Monitoring of Temperature. When you wish to monitor the actual temperature of a distillation, a thermometer must be placed inside the apparatus. Figures 14.7 and 14.8 show distillation assemblies that use an internal thermometer. The apparatus in Figure 14.7A represents the simplest possible distillation assembly. It does not use a condenser, and the thermometer is suspended from a clamp. It is possible to add either an air or a water condenser to this basic assembly (Figure 14.7B) and maintain internal monitoring of the temperature.

In the arrangement shown in Figure 14.8, a thermometer adapter is used. A thermometer adapter (Figure 14.9A) provides a convenient way of holding a thermometer in place. The Claisen head is used to provide an opening to the atmosphere, thereby avoiding a sealed system. With the Claisen head, a drying tube may be used to protect the system from atmospheric moisture.

If protection from atmospheric moisture is not required, the multipurpose adapter may be used. The multipurpose adapter (Figure 14.9B) replaces both the thermometer adapter and the Claisen head. With this adapter, the necessary opening to the atmosphere is provided by the side arm. The threaded joint holds the thermometer in place.

Carefully notice the position of the thermometer in Figures 14.7 and 14.8. The bulb of the thermometer must be placed in the stem of the Hickman head, *just below*

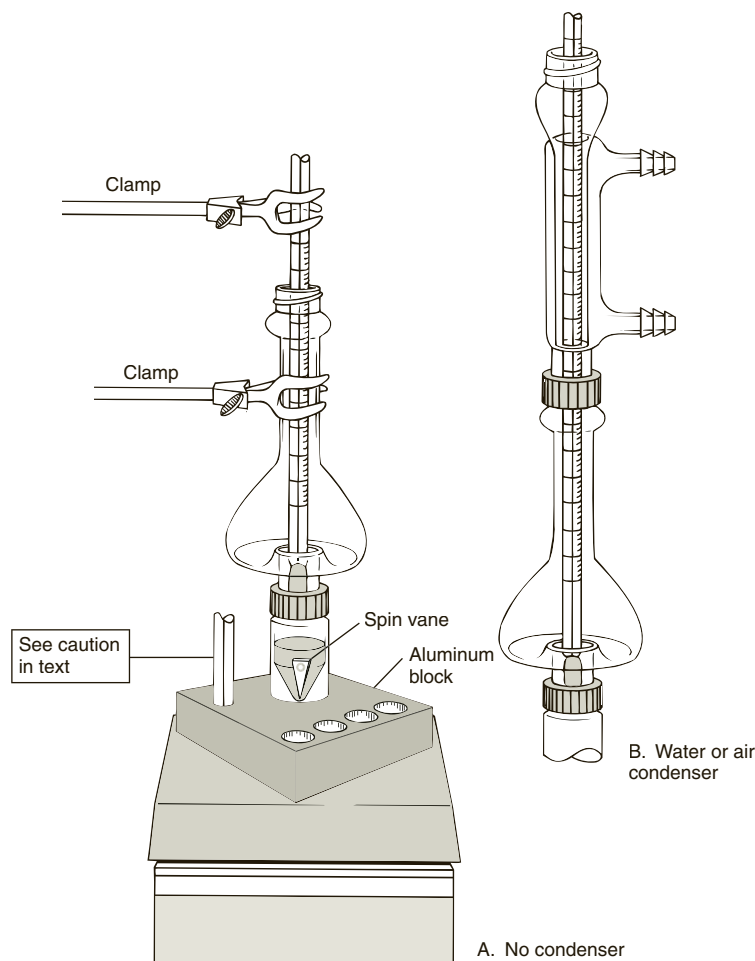


Figure 14.7
Basic microscale distillation (internal monitoring of temperature). Do not use a mercury thermometer.

the well, or it will not read the temperature correctly. The distillation temperature can be monitored most accurately by using a partial immersion mercury thermometer (see Technique 13, Section 13.3).

NOTE: It is good practice to monitor the temperature internally whenever possible.

Boiling Stones or Stirring. A boiling stone should be used during distillation in order to prevent bumping. As an alternative, the liquid being distilled may be rapidly stirred. A triangular spin vane of the correct size should be used when distilling from a conical vial, whereas a stirring bar should be used when distilling from a round-bottom flask.

Size of Distillation Flask. As a rule, the distillation flask or vial should not be filled to more than two thirds of its total capacity. This allows room for boiling and stirring action, and it prevents contamination of the distillate by bumping. A flask that is too large should also be avoided. With too large a flask, the **holdup** is excessive; the holdup is the amount of material that cannot distill because some vapor must fill the empty flask.

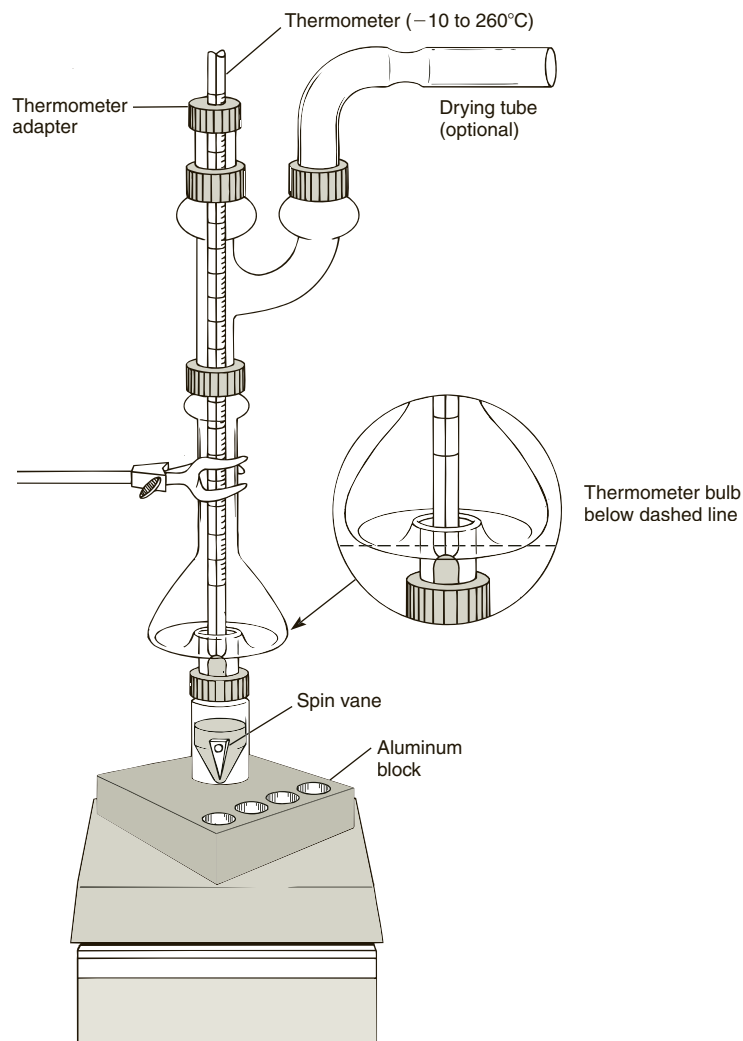


Figure 14.8
Basic microscale distillation using thermometer adapter
(internal monitoring of temperature).

Assembling the Apparatus. You should not grease the joints when assembling the apparatus. Ungreased joints seal well enough to allow you to perform a simple distillation. Stopcock grease can introduce a serious contaminant into your product.

Rate and Degree of Heating. Take care not to distill too quickly. If you vaporize liquid at a rate faster than it can be recondensed, some of your product may be lost by evaporation. On the other hand, you should not distill too slowly. This may also lead to loss of product because there is a longer period during which vapors can escape. Carefully examine your apparatus during distillation to monitor the position of either a reflux ring or a wet appearance on the surface of the glass. Either of these indicates the place at which condensation is occurring. The position at which condensation occurs should be well inside the Hickman head. Be sure that liquid is collecting in the well. If all the surfaces are shiny (wet) and there is no distillate, you are losing material.

NOTE: A slower rate of heating also helps to avoid bumping.

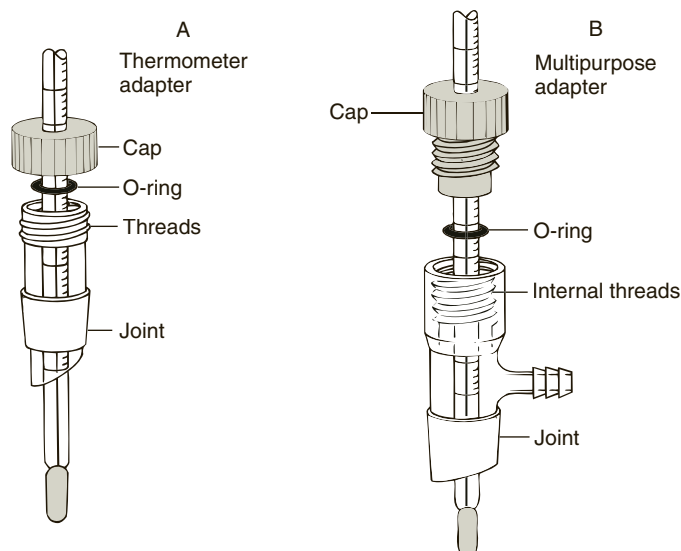


Figure 14.9
Two adapters.

If you are using a sand bath, material may be lost because the hot sand bath radiates too much heat upward and warms the Hickman still. If you believe this to be the case, it can often be remedied by placing a small square of aluminum foil over the top of the sand bath. Make a tear from one edge to the center of the foil to wrap it around the apparatus.

14.4 Semi-microscale and Macroscale Equipment

When you wish to distill quantities of liquid that are larger than 2–3 mL, different equipment is required. Most manufacturers of microscale equipment make two pieces of conventional distillation equipment sized to work with the $\text{F } 14/10$ microscale kit components. These two pieces, the **distillation head** and the **bent vacuum adapter**, are not provided in student microscale kits but must be purchased separately. Figure 14.10 shows a semimicroscale assembly using these components. Note that the bulb of the thermometer must be placed *below the side arm* if it is to be bathed in vapor and give a correct temperature reading. This apparatus assumes that a condenser is not necessary; however, you could easily insert one between the distilling head and the bent vacuum adapter. This insertion would produce a completely traditional distillation apparatus but would use microscale equipment. A distillation apparatus constructed from a “macroscale” organic laboratory kit is shown in Figure 14.11. This type of equipment is being used today in organic laboratories that have not converted to microscale. Electrically regulated **heating mantles** are often used with this equipment.

14.5 Using a Digital Thermometer

Rather than using a glass thermometer to monitor the temperature during a distillation, one can use a digital thermometer with a stainless steel temperature probe. The use of a digital thermometer for a semimicroscale distillation is illustrated in Figure 14.12. Digital thermometers may, of course, also be used in other styles of apparatus such as a standard distillation with a condenser or a microscale distillation. Technique 13, Section 13.4 discusses the use of digital thermometers, or the LabQuest and LabPro systems, in more detail.

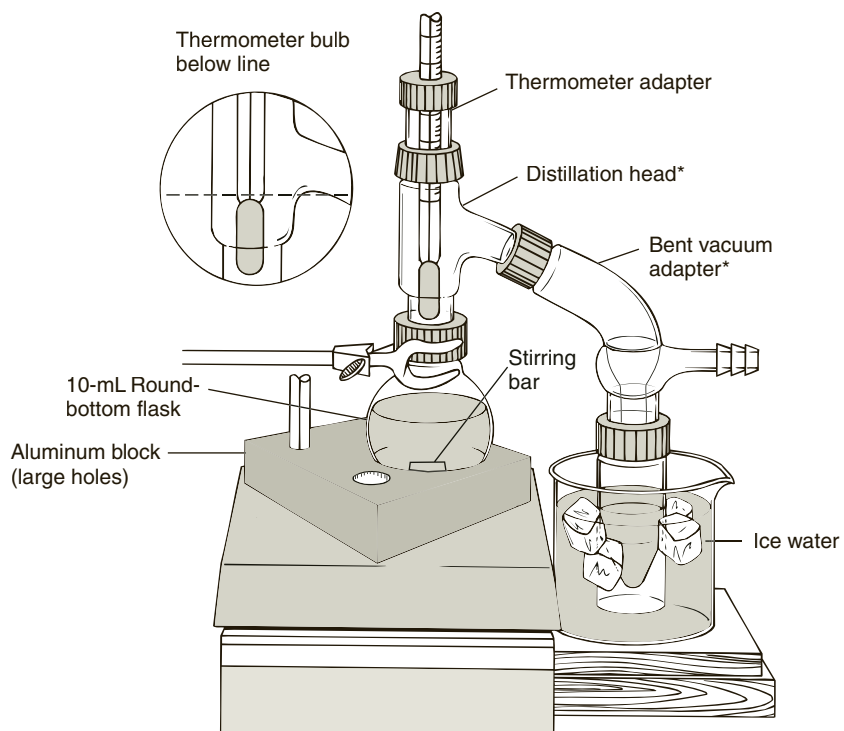


Figure 14.10
Semimicroscale distillation (*requires special pieces).

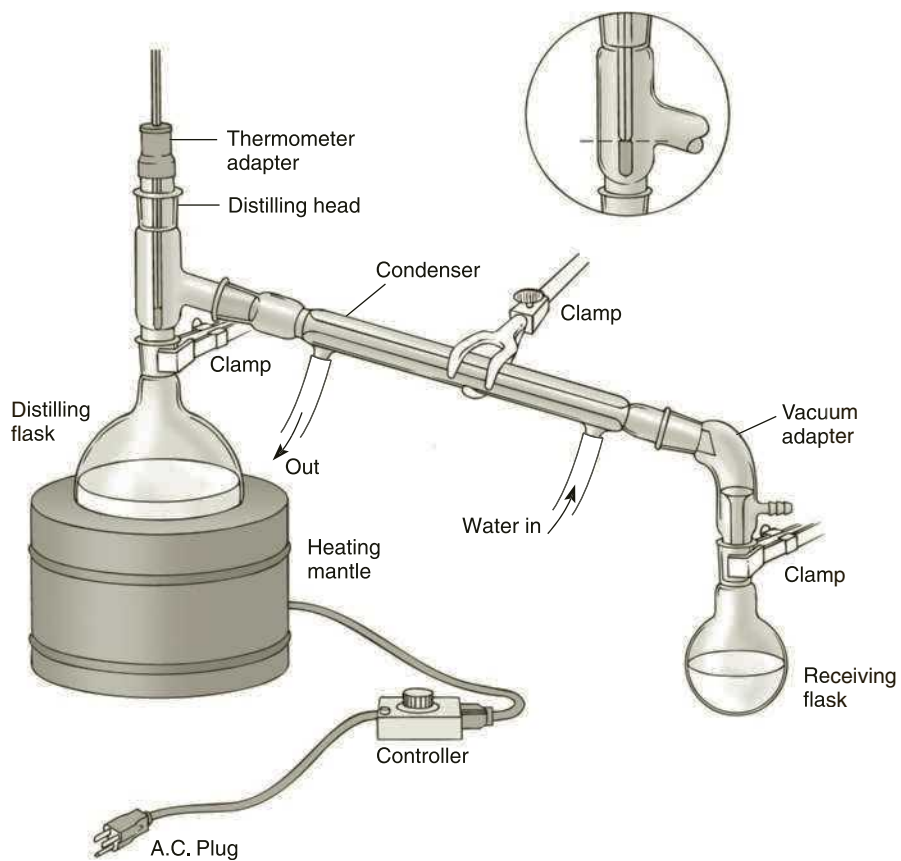


Figure 14.11
Distillation with the standard macroscale organic lab kit.

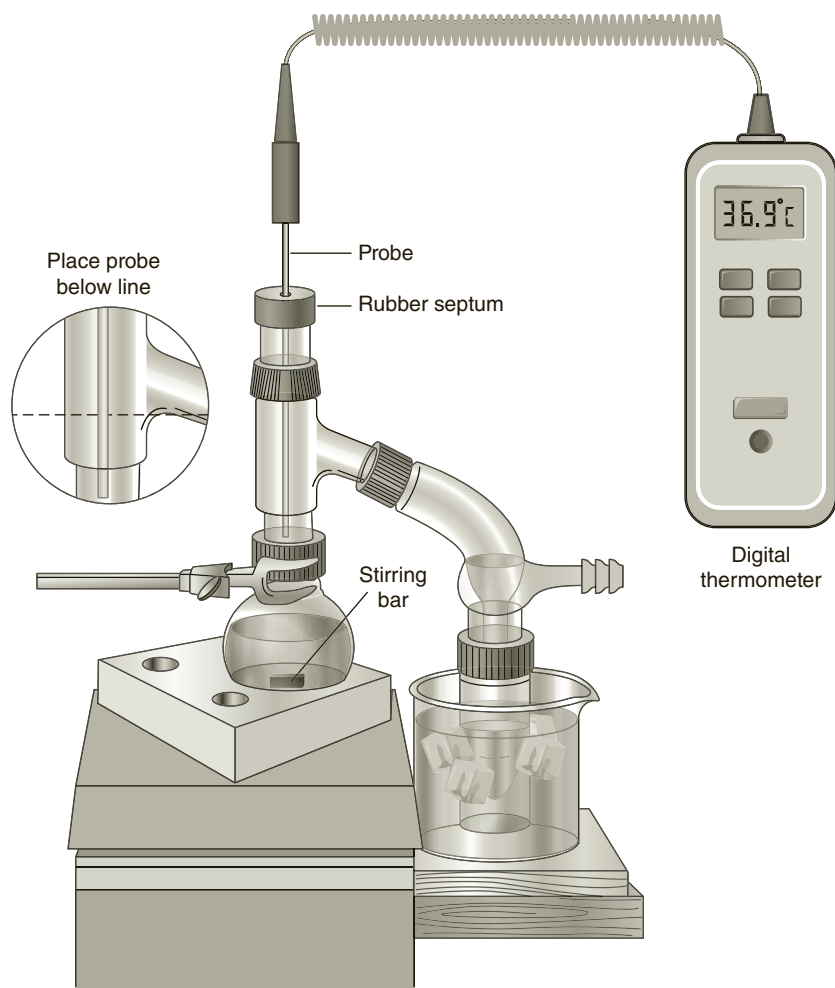


Figure 14.12
Distillation using a digital thermometer.

PROBLEMS

- Using Figure 14.3, answer the following questions.
 - What is the molar composition of the vapor in equilibrium with a boiling liquid that has a composition of 60% A and 40% B?
 - A sample of vapor has the composition 50% A and 50% B. What is the composition of the boiling liquid that produced this vapor?
- Use an apparatus similar to that shown in Figure 14.10 and assume that the round-bottom flask holds 10 mL and that the Claisen head has an internal volume of about 2 mL in the vertical section. At the end of a distillation, vapor would fill this volume, but it could not be forced through the system. No liquid would remain in the distillation flask. Assuming this **holdup volume** of 12 mL, use the ideal gas law and assume a boiling point of 100°C (760 mm Hg) to calculate the number of microliters of liquid ($d = 0.9 \text{ g/mL}$, $MW = 200$) that would recondense into the distillation flask on cooling.
- Explain the significance of a horizontal line connecting a point on the lower curve with a point on the upper curve (such as line XY) in Figure 14.3.

4. Using Figure 14.3, determine the boiling point of a liquid having a molar composition of 50% A and 50% B.
5. What is the approximate difference between the temperature of a boiling liquid in a conical vial and the temperature read on an *external* thermometer when both are placed on an aluminum block?
6. Where should the thermometer bulb be located for internal monitoring in
 - a. a distillation apparatus using a Hickman head?
 - b. a large-scale distillation using a Claisen head with a water condenser placed beyond it?
7. Under what conditions can a good separation be achieved with a simple distillation?

Fractional Distillation, Azeotropes

Simple distillation, described in Technique 14, works well for most routine separation and purification procedures for organic compounds. When boiling-point differences of components to be separated are not large, however, **fractional distillation** must be used to achieve a good separation.

PART A. FRACTIONAL DISTILLATION

15.1 Differences between Simple and Fractional Distillation

When an ideal solution of two liquids, such as benzene (bp 80°C) and toluene (bp 110°C), is distilled by simple distillation, the first vapor produced will be enriched in the lower-boiling component (benzene). However, when that initial vapor is condensed and analyzed, the distillate will not be pure benzene. The boiling-point difference of benzene and toluene (30°C) is too small to achieve a complete separation by simple distillation. Following the principles outlined in Technique 14, Section 14.2, and using the vapor–liquid composition curve given in Figure 15.1, you can see what would happen if you started with an equimolar mixture of benzene and toluene.

Following the dashed lines shows that an equimolar mixture (50 mole % benzene) would begin to boil at about 91°C and, far from being 100% benzene, the distillate would contain about 74 mole % benzene and 26 mole % toluene. As the distillation continued, the composition of the undistilled liquid would move in the direction of A' (there would be increased toluene due to removal of more benzene than toluene), and the corresponding vapor would contain a progressively smaller amount of benzene. In effect, the temperature of the distillation would continue to increase throughout the distillation (as in Technique 14, Figure 14.2B), and it would

be impossible to obtain any fraction that consisted of pure benzene.

Suppose, however, that we are able to collect a small quantity of the first distillate that was 74 mole % benzene, and to redistill it. Using Figure 15.1, we can see that this liquid would begin to boil at about 84°C and would give an initial distillate containing 90 mole % benzene. If we were experimentally able to continue taking small fractions at the beginning of each distillation and redistill them, we would eventually reach a liquid with a composition of nearly 100 mole % benzene. However, because we only took a small amount of material at the beginning of each distillation, we would have lost most of the material we started with. To recapture a reasonable amount of benzene, we would have to process each of the fractions left behind in

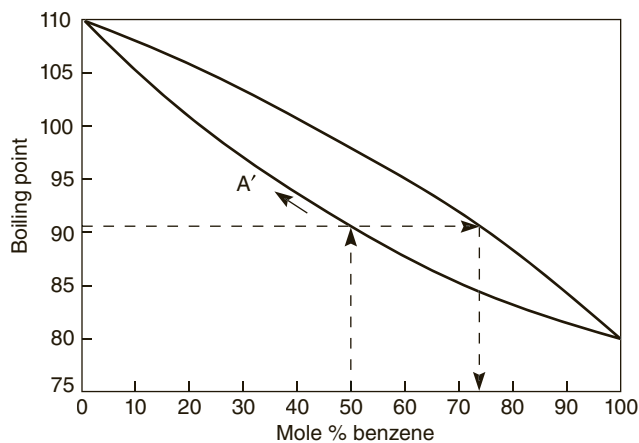


Figure 15.1
The vapor–liquid composition curve for mixtures of benzene and toluene.

the same way as our early fractions. As each of them was partially distilled, the material advanced would become progressively richer in benzene, whereas that left behind would become progressively richer in toluene. It would require thousands (maybe millions) of such microdistillations to separate benzene from toluene.

Obviously, the procedure just described would be very tedious; fortunately, it need not be performed in usual laboratory practice. **Fractional distillation** accomplishes the same result. You simply have to use a column inserted between the distillation flask and the receiver (Hickman head), as shown in Figure 15.2. This **fractionating column** is filled, or **packed**, with a suitable material such as a stainless steel sponge. This packing allows a mixture of benzene and toluene to be subjected continuously to many vaporization–condensation cycles as the material moves up the column. With each cycle within the column, the composition of the vapor is progressively enriched in the lower-boiling component (benzene). Nearly pure benzene (bp 80°C) finally emerges from the top of the column, condenses, and passes into the receiving head or flask. This process continues until all the benzene is removed. The distillation must be carried out slowly to ensure that numerous

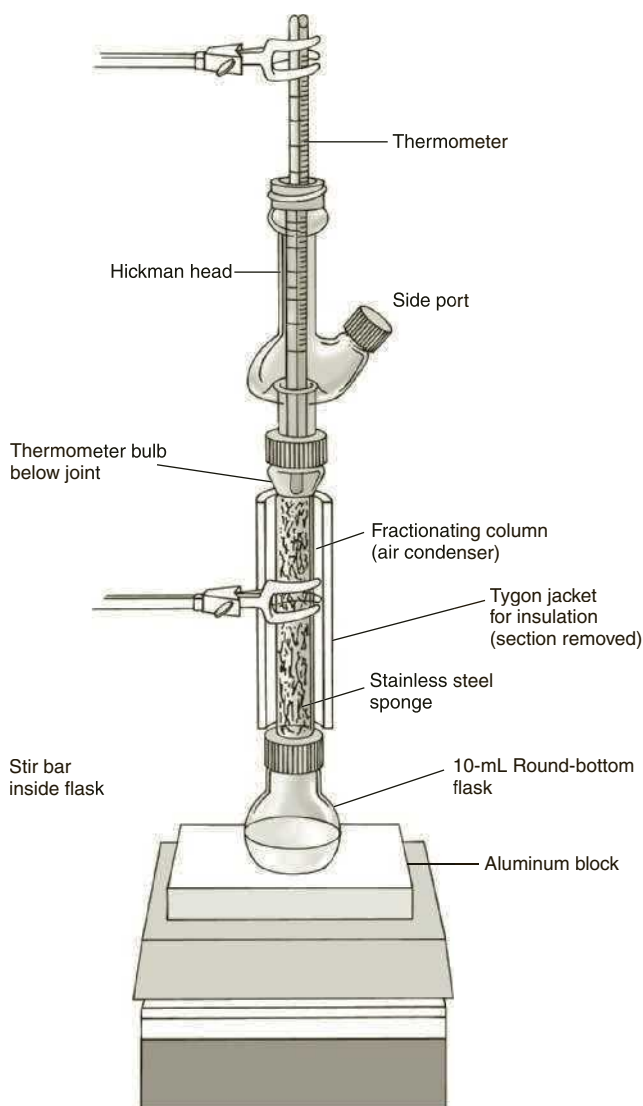
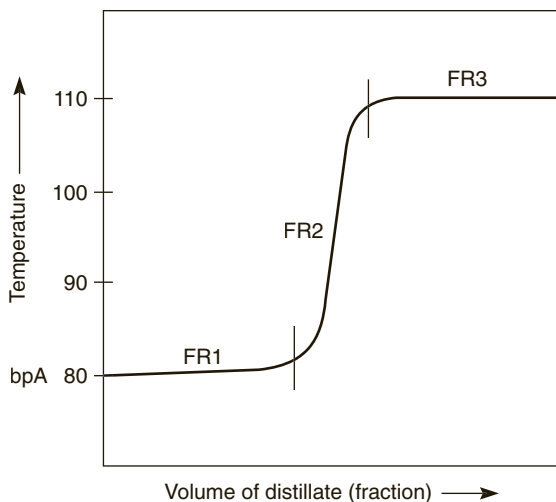
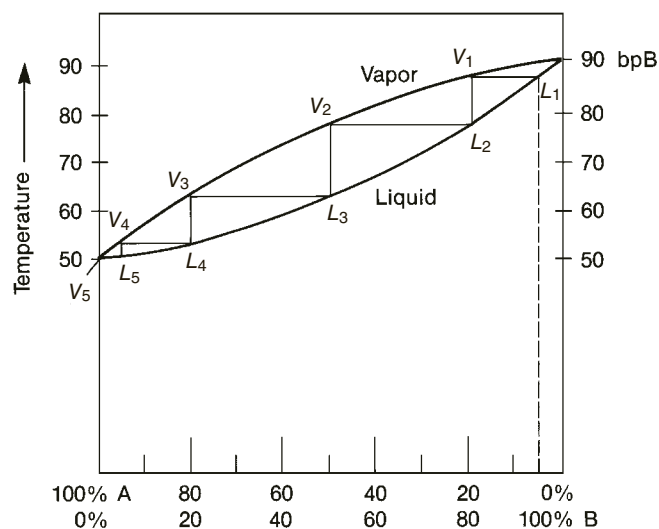


Figure 15.2
Microscale apparatus for fractional distillation.

**Figure 15.3**

Temperature–distillate plot for fractional distillation of a benzene–toluene mixture.

**Figure 15.4**

Phase diagram for a fractional distillation of an ideal two-component system.

vaporization–condensation cycles occur. When nearly all the benzene has been removed, the temperature begins to rise, and a small amount of a second fraction, which contains some benzene and toluene, may be collected. When the temperature reaches 110°C, the boiling point of pure toluene, the vapor is condensed and collected as the third fraction. A plot of boiling point versus volume of condensate (distillate) would resemble Figure 15.3. This separation would be much better than that achieved by simple distillation (Figure 15.1).

15.2 Vapor–Liquid Composition Diagrams

A vapor–liquid composition–phase diagram like the one in Figure 15.4 can be used to explain the operation of a fractionating column with an **ideal solution** of two liquids, A and B. An ideal solution is one in which the two liquids are chemically similar, miscible (mutually soluble) in all proportions, and do not interact. Ideal solutions obey **Raoult's Law**. Raoult's Law is explained in detail in Section 15.3.

The phase diagram relates the compositions of the boiling liquid (lower curve) and its vapor (upper curve) as a function of temperature. Any horizontal line drawn across the diagram (a constant-temperature line) intersects the diagram in two places. These intersections relate the vapor composition to the composition of the boiling liquid that produces that vapor. By convention, composition is expressed either in **mole fraction** or in **mole percentage**. The mole fraction is defined as follows:

$$\text{Mole fraction A} = N_A = \frac{\text{moles A}}{\text{moles A} + \text{moles B}}$$

$$\text{Mole fraction B} = N_B = \frac{\text{moles B}}{\text{moles A} + \text{moles B}}$$

$$N_A + N_B = 1$$

$$\text{Mole percentage A} = N_A \times 100$$

$$\text{Mole percentage B} = N_B \times 100$$

The horizontal and vertical lines shown in Figure 15.4 represent the processes that occur during a fractional distillation. Each of the **horizontal lines** (L_1V_1 , L_2V_2 , etc.) represents the **vaporization** step of a given vaporization–condensation cycle

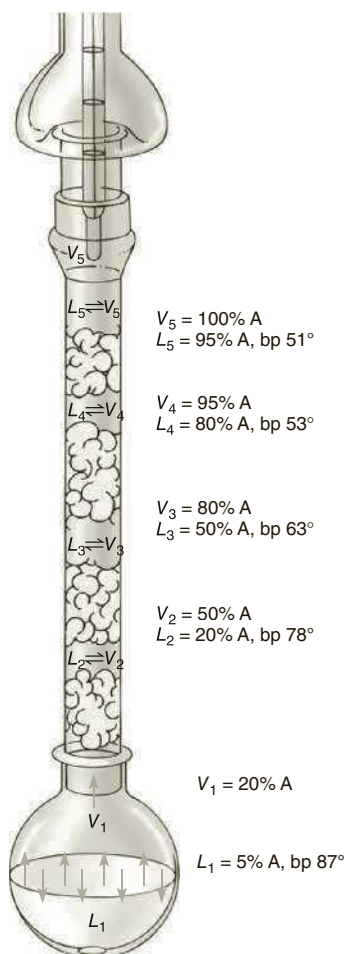


Figure 15.5
Vaporization–condensation
in a fractionation column.

and represents the composition of the vapor in equilibrium with liquid at a given temperature. For example, at 63°C a liquid with a composition of 50% A (L_3 on the diagram) would yield vapor of composition 80% A (V_3 on diagram) at equilibrium. The vapor is richer in the lower-boiling component A than the original liquid was.

Each of the **vertical lines** (V_1L_2 , V_2L_3 , etc.) represents the **condensation** step of a given vaporization–condensation cycle. The composition does not change as the temperature drops on condensation. The vapor at V_3 , for example, condenses to give a liquid (L_4 on the diagram) of composition 80% A with a drop in temperature from 63° to 53°C.

In the example shown in Figure 15.4, pure A boils at 50°C and pure B boils at 90°C. These two boiling points are represented at the left- and right-hand edges of the diagram, respectively. Now consider a solution that contains only 5% of A but 95% of B. (Remember that these are *mole* percentages.) This solution is heated (following the dashed line) until it is observed to boil at L_1 (87°C). The resulting vapor has composition V_1 (20% A, 80% B). The vapor is richer in A than the original liquid, but it is by no means pure A. In a simple distillation apparatus, this vapor would be condensed and passed into the receiver in a very impure state. However, with a fractionating column in place, the vapor is condensed in the **column** to give liquid L_2 (20% A, 80% B). Liquid L_2 is immediately revaporized (bp 78°C) to give a vapor of composition V_2 (50% A, 50% B), which is condensed to give liquid L_3 . Liquid L_3 is revaporized (bp 63°C) to give vapor of composition V_3 (80% A, 20% B), which is condensed to give liquid L_4 . Liquid L_4 is revaporized (bp 53°C) to give vapor of composition V_4 (95% A, 5% B). This process continues to V_5 , which condenses to give nearly pure liquid A. The fractionating process follows the stepped lines in the figure downward and to the left.

As this process continues, all of liquid A is removed from the distillation flask or vial, leaving nearly pure B behind. If the temperature is raised, liquid B may be distilled as a nearly pure fraction. Fractional distillation will have achieved a separation of A and B, a separation that would have been nearly impossible with simple distillation. Notice that the boiling point of the liquid becomes lower each time it vaporizes. Because the temperature at the bottom of a column is normally higher than the temperature at the top, successive vaporizations occur higher and higher in the column as the composition of the distillate approaches that of pure A. This process is illustrated in Figure 15.5, where the composition of the liquids, their boiling points, and the composition of the vapors present are shown alongside the fractionating column.

15.3 Raoult's Law

Two liquids (A and B) that are miscible and that do not interact form an **ideal solution** and follow Raoult's Law. The law states that the partial vapor pressure of component A in the solution (P_A) equals the vapor pressure of pure A (P_A°) times its mole fraction (N_A) (Equation 1). A similar expression can be written for component B (Equation 2). The mole fractions N_A and N_B were defined in Section 15.2.

$$\text{Partial vapor pressure of A in solution} = P_A = (P_A^\circ)(N_A) \quad [1]$$

$$\text{Partial vapor pressure of B in solution} = P_B = (P_B^\circ)(N_B) \quad [2]$$

P_A° is the vapor pressure of pure A, independent of B. P_B° is the vapor pressure of B, independent of A. In a mixture of A and B, the partial vapor pressures are added to give the total vapor pressure above the solution (Equation 3). When the total pressure (sum of the partial pressures) equals the applied pressure, the solution boils.

$$P_{\text{total}} = P_A + P_B = P_A^\circ N_A + P_B^\circ N_B \quad [3]$$

The composition of A and B in the vapor produced is given by Equations 4 and 5.

$$N_A(\text{vapor}) = \frac{P_A}{P_{\text{total}}} \quad [4]$$

$$N_B(\text{vapor}) = \frac{P_B}{P_{\text{total}}} \quad [5]$$

Several problems involving applications of Raoult's Law are illustrated in Figure 15.6. Note, particularly in the result from Equation 4, that the vapor is richer ($N_A = 0.67$) in the lower-boiling (higher vapor pressure) component A than it was before vaporization ($N_A = 0.50$). This proves mathematically what was described in Section 15.2.

The consequences of Raoult's Law for distillations are shown schematically in Figure 15.7. In Part A, the boiling points are identical (vapor pressures the same), and no separation is attained regardless of how the distillation is conducted. In Part B, a fractional distillation is required, whereas in Part C a simple distillation provides an adequate separation.

When a solid B (rather than another liquid) is dissolved in a liquid A, the boiling point is increased. In this extreme case, the vapor pressure of B is negligible, and the vapor will be pure A no matter how much solid B is added. Consider a solution of salt in water.

$$\begin{aligned} P_{\text{total}} &= P_{\text{water}}^\circ N_{\text{water}} + P_{\text{salt}}^\circ N_{\text{salt}} \\ P_{\text{salt}}^\circ &= 0 \\ P_{\text{total}} &= P_{\text{water}}^\circ N_{\text{water}} \end{aligned}$$

A solution whose mole fraction of water is 0.7 will not boil at 100°C, because $P_{\text{total}} = (760)(0.7) = 532$ mm Hg and is less than atmospheric pressure. If the solution is heated to 110°C, it will boil because $P_{\text{total}} = (1085)(0.7) = 760$ mm Hg. Although the solution must be heated to 110°C to boil it, the vapor is pure water and has a boiling-point temperature of 100°C. (The vapor pressure of water at 110°C can be looked up in a handbook; it is 1085 mm Hg.)

Consider a solution at 100°C where $N_A = 0.5$ and $N_B = 0.5$.

1. What is the partial vapor pressure of A in the solution if the vapor pressure of pure A at 100°C is 1020 mm Hg?

Answer: $P_A = P_A^\circ N_A = (1020)(0.5) = 510$ mm Hg

2. What is the partial vapor pressure of B in the solution if the vapor pressure of pure B at 100°C is 500 mm Hg?

Answer: $P_B = P_B^\circ N_B = (500)(0.5) = 250$ mm Hg

3. Would the solution boil at 100°C if the applied pressure were 760 mm Hg?

Answer: Yes. $P_{\text{total}} = P_A + P_B = (510 + 250) = 760$ mm Hg


4. What is the composition of the vapor at the boiling point?

Answer: The boiling point is 100°C

$$N_A(\text{vapor}) = \frac{P_A}{P_{\text{total}}} = 510/760 = 0.67$$

$$N_B(\text{vapor}) = \frac{P_B}{P_{\text{total}}} = 250/760 = 0.33$$

Figure 15.6
Sample calculations with Raoult's Law.

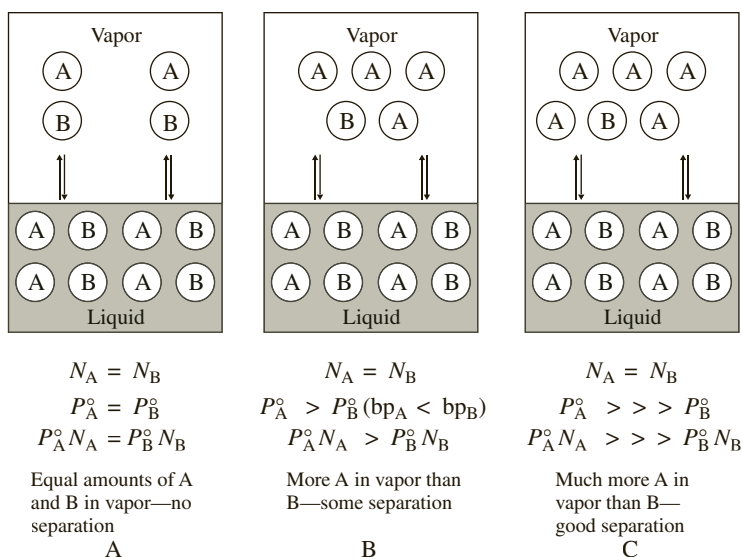**Figure 15.7**

Consequences of Raoult's Law. (A) Boiling points (vapor pressures) are identical—no separation. (B) Boiling point somewhat lower for A than for B—requires fractional distillation. (C) Boiling point much lower for A than for B—simple distillation will suffice.

15.4 Column Efficiency

A common measure of the efficiency of a column is given by its number of **theoretical plates**. The number of theoretical plates in a column is related to the number of vaporization–condensation cycles that occur as a liquid mixture travels through it. Using the example mixture in Figure 15.4, if the first distillate (condensed vapor) had the composition at L_2 when starting with liquid of composition L_1 , the column would be said to have *one theoretical plate*. This would correspond to a simple distillation, or one vaporization–condensation cycle. A column would have two theoretical plates if the first distillate had the composition at L_3 . The two-theoretical-plate column essentially carries out “two simple distillations.” According to Figure 15.4, *five theoretical plates* would be required to separate the mixture that started with composition L_1 . Notice that this corresponds to the number of “steps” that need to be drawn in the figure to arrive at a composition of 100% A.

Most columns do not allow distillation in discrete steps, as indicated in Figure 15.4. Instead, the process is *continuous*, allowing the vapors to be continuously in contact with liquid of changing composition as they pass through the column. Any material can be used to pack the column as long as it can be wetted by the liquid and as long as it does not pack so tightly that vapor cannot pass.


The approximate relationship between the number of theoretical plates needed to separate an ideal two-component mixture and the difference in boiling points is given in Table 15.1. Notice that more theoretical plates are required as the boiling-point differences between the components decrease. For instance, a mixture of A (bp 130°C) and B (bp 166°C) with a boiling-point difference of 36°C would be expected to require a column with a minimum of five theoretical plates.

15.5 Types of Fractionating Columns and Packings

Several types of fractionating columns are shown in Figure 15.8. The Vigreux column, shown in Part A, has indentations that incline downward at angles of 45 degrees and are in pairs on opposite sides of the column. The projections into the column provide increased possibilities for condensation and for the vapor to

Table 15.1 Theoretical plates required to separate mixtures, based on boiling-point differences of components

<i>Boiling-Point Difference</i>	<i>Number of Theoretical Plates</i>
108	1
72	2
54	3
43	4
36	5
20	10
10	20
7	30
4	50
2	100

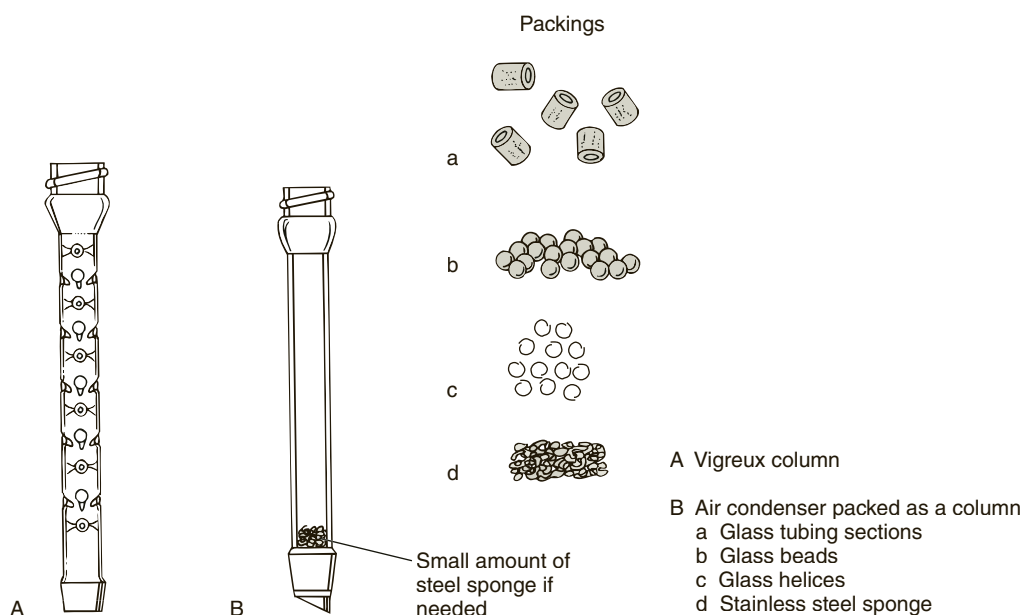**Figure 15.8**
Columns for fractional distillation.

equilibrate with the liquid. Vigreux columns are popular in cases where only a small number of theoretical plates are required. They are not very efficient (a 20-cm column might have only 2.5 theoretical plates), but they allow for rapid distillation and have a small **holdup** (the amount of liquid retained by the column). A column packed with a stainless steel sponge is a more effective fractionating column than a Vigreux column, but not by a large margin. Glass beads, or glass helices, can also be used as a packing material, and they have a slightly greater efficiency. The air condenser or the water condenser can be used as an improvised column if an actual fractionating column is unavailable. If a condenser is packed with glass beads, glass helices, or sections of glass tubing, the packing must be held in place by inserting a small plug of stainless steel sponge into the bottom of the condenser.

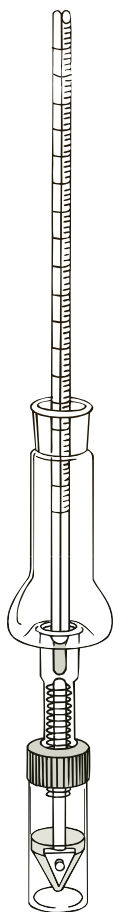


Figure 15.10
A commercially available microscale spinning-band column.

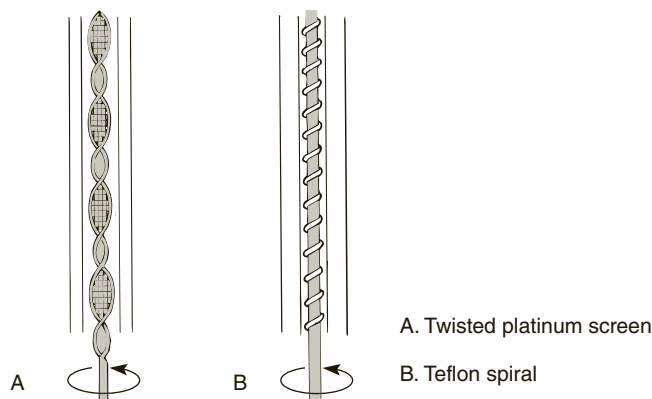


Figure 15.9
Bands for spinning-band columns.

The most effective type of column is the **spinning-band column**. In the most elegant form of this device, a tightly fitting, twisted platinum screen or a Teflon rod with helical threads is rotated rapidly inside the bore of the column (Figure 15.9). A spinning-band column that is available for microscale work is shown in Figure 15.10. This spinning-band column has a band about 2–3 cm in length and provides 4–5 theoretical plates. It can separate 1–2 mL of a mixture with a 30°C boiling-point difference. Larger research models of this spinning-band column can provide as many as 20 or 30 theoretical plates and can separate mixtures with a boiling-point difference of as little as 5–10°C.

Manufacturers of fractionating columns often offer them in a variety of lengths. Because the efficiency of a column is a function of its length, longer columns have more theoretical plates than shorter ones. It is common to express efficiency of a column in a unit called **HETP**, the **H**eight of a column that is **E**quivalent to one **T**heoretical **P**late. HETP is usually expressed in units of cm/plate. When the height of the column (in centimeters) is divided by this value, the total number of theoretical plates is specified.

15.6 Fractional Distillation: Methods and Practice

When a fractional distillation is performed, the column should be clamped in a vertical position. The distillation should be conducted as slowly as possible, but the rate of distillation should be steady enough to produce a constant temperature reading at the thermometer.

Many fractionating columns must be insulated so that temperature equilibrium is maintained at all times. Additional insulation will not be required for columns that have an evacuated outer jacket, but those that do not can benefit from being wrapped in insulation.

A microscale air condenser can be converted to a column by packing it with a piece of stainless steel sponge. The simplest form of insulation is Tygon tubing that has been split lengthwise. Select a piece with an inner diameter that just matches or is slightly smaller than the diameter of the fractionating column so that it will fit snugly.

CAUTION



Cut the tubing to the correct length and then slit it with a sharp scissors. Do not use a razor blade or knife. Tygon tubing is difficult to cut; it is a nonslip substance and will “grab” even a single-edged razor blade in a way that can give you a nasty cut. See Experiment 6, page 52, for complete instructions.

Glass wool and aluminum foil (shiny side in) are often used for insulation. You can wrap the column with glass wool and then use a wrapping of the aluminum foil to keep it in place. An especially effective method is to make an insulation blanket by placing a layer of glass wool or cotton between two rectangles of aluminum foil, placed shiny side in. The sandwich is bound together with duct tape. This blanket, which is reusable, can be wrapped around the column and held in place with twist ties or tape.

The **reflux ratio** is defined as the ratio of the number of drops of distillate that return to the distillation flask compared to the number of drops of distillate collected. In an efficient column, the reflux ratio should equal or exceed the number of theoretical plates. A high reflux ratio ensures that the column will achieve temperature equilibrium and achieve its maximum efficiency. This ratio is not easy to determine; in fact, it is impossible to determine when using a Hickman head, and it should not concern a beginning student. In some cases, the **throughput**, or **rate of takeoff**, of a column may be specified. This is expressed as the number of milliliters of distillate that can be collected per unit of time, usually as mL/min.

Microscale Apparatus. The apparatus shown in Figure 15.2 is the one you are most likely to use in the microscale laboratory. If your laboratory is one of the better-equipped ones, you may have access to spinning-band columns like those shown in Figure 15.10. The distillation temperature can be monitored most accurately by using a partial immersion mercury thermometer (see Technique 13, Section 13.3).

Macroscale Apparatus. Figure 15.11 illustrates a fractional distillation assembly that can be used for larger-scale distillations. It has a glass-jacketed column that is packed with a stainless steel sponge. This apparatus would be common in situations where quantities of liquid in excess of 10 mL were to be distilled.

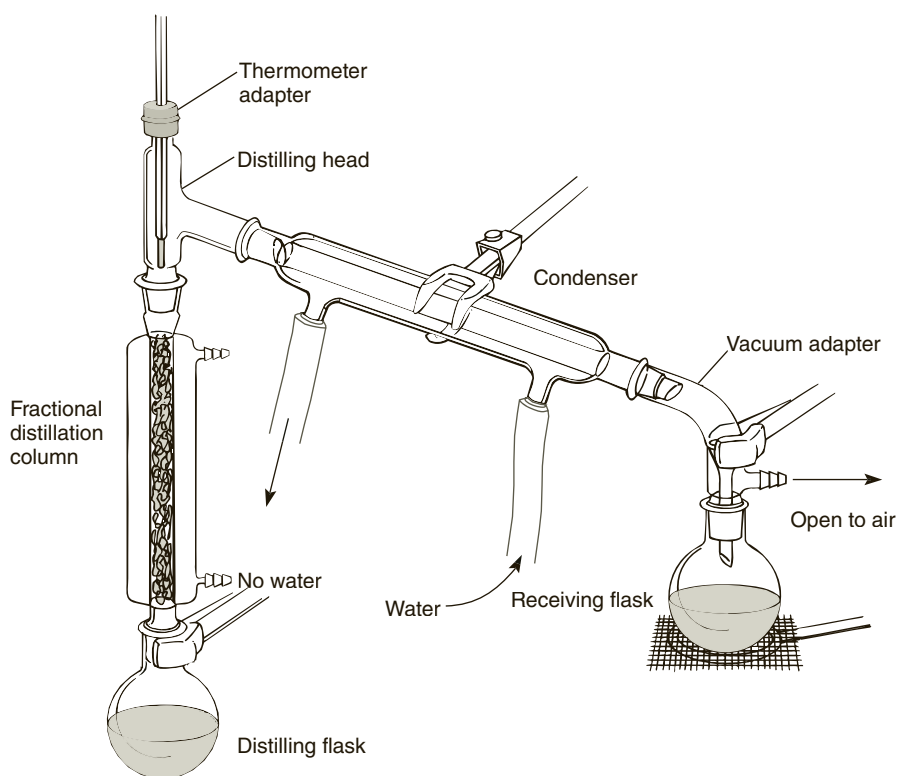


Figure 15.11
Large-scale fractional distillation apparatus.

PART B. AZEOTROPES

15.7 Nonideal Solutions: Azeotropes

Some mixtures of liquids, because of attractions or repulsions between the molecules, do not behave ideally; they do not follow Raoult's Law. There are two types of vapor–liquid composition diagrams that result from this nonideal behavior: **minimum-boiling-point** and **maximum-boiling-point** diagrams. The minimum or maximum points in these diagrams correspond to a constant-boiling mixture called an **azeotrope**. An azeotrope is a mixture with a fixed composition that cannot be altered by either simple or fractional distillation. An azeotrope behaves as if it were a pure compound, and it distills from the beginning to the end of its distillation at a constant temperature, giving a distillate of constant (azeotropic) composition. The vapor in equilibrium with an azeotropic liquid has the same composition as the azeotrope. Because of this, an azeotrope is represented as a *point* on a vapor–liquid composition diagram.

A. Minimum-Boiling-Point Diagrams

A minimum-boiling-point azeotrope results from a slight incompatibility (repulsion) between the liquids being mixed. This incompatibility leads to a higher-than-expected combined vapor pressure from the solution. This higher combined vapor pressure brings about a lower boiling point for the mixture than is observed for the pure components. The most common two component mixture that gives a minimum-boiling-point azeotrope is the ethanol–water system shown in Figure 15.12. The azeotrope at V_3 has a composition of 96% ethanol–4% water and a boiling point of 78.1°C . This boiling point is not much lower than that of pure ethanol (78.3°C), but it means that it is impossible to obtain pure ethanol from the distillation of any ethanol–water mixture that contains more than 4% water. Even with the best fractionating column, you cannot obtain 100% ethanol. The remaining 4% of water can be removed by adding benzene and removing a different azeotrope, the ternary benzene–water–ethanol azeotrope (bp 65°C). Once the water is removed, the excess benzene is removed as an ethanol–benzene azeotrope (bp 68°C). The resulting material is free of water and is called “absolute” ethanol.

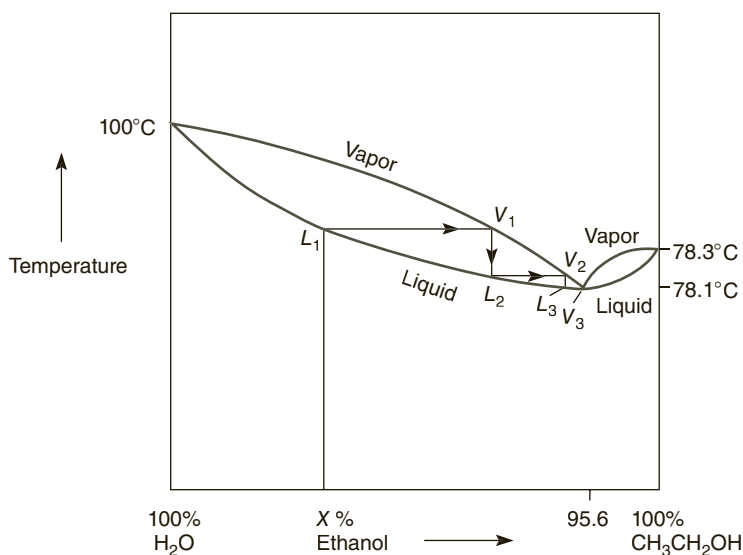


Figure 15.12
Ethanol–water minimum-boiling-point phase diagram.

The fractional distillation of an ethanol–water mixture of composition X can be described as follows. The mixture is heated (follow line XL_1) until it is observed to boil at L_1 . The resulting vapor at V_1 will be richer in the lower-boiling component, ethanol, than the original mixture.¹ The condensate at L_2 is vaporized to give V_2 . The process continues, following the lines to the right, until the azeotrope is obtained at V_3 . The liquid that distills is not pure ethanol, but it has the azeotropic composition of 96% ethanol and 4% water, and it distills at 78.1°C. The azeotrope, which is richer in ethanol than the original mixture, continues to distill. As it distills, the percentage of water left behind in the distillation flask continues to increase. When all the ethanol has been distilled (as the azeotrope), pure water remains behind in the distillation flask, and it distills at 100°C.

If the azeotrope obtained by the preceding procedure is redistilled, it distills from the beginning to the end of the distillation at a constant temperature of 78.1°C as if it were a pure substance. There is no change in the composition of the vapor during the distillation.

Some common minimum-boiling azeotropes are given in Table 15.2. Numerous other azeotropes are formed in two- and three-component systems; such azeotropes are common. Water forms azeotropes with many substances; therefore, water must be carefully removed with **drying agents** whenever possible before compounds are distilled. Extensive azeotropic data are available in references such as the *Handbook of Chemistry and Physics*.²

Note that azeotropic combinations with water should be miscible or at least have a high degree of mutual solubility. If the water does not dissolve in the other component, or has low solubility, then the process becomes a steam distillation (see Technique 18).

B. Maximum-Boiling-Point Diagrams

A maximum-boiling-point azeotrope results from a slight attraction between the component molecules. This attraction leads to lower combined vapor pressure than

Table 15.2 Common minimum-boiling-point azeotropes

Azeotrope	Composition (Weight Percentage)	Boiling Point (°C)
Ethanol–water	95.6% C ₂ H ₅ OH, 4.4% H ₂ O	78.17
Benzene–water	91.1% C ₆ H ₆ , 8.9% H ₂ O	69.4
Benzene–water–ethanol	74.1% C ₆ H ₆ , 7.4% H ₂ O, 18.5% C ₂ H ₅ OH	64.9
Methanol–carbon tetrachloride	20.6% CH ₃ OH, 79.4% CCl ₄	55.7
Ethanol–benzene	32.4% C ₂ H ₅ OH, 67.6% C ₆ H ₆	67.8
Methanol–toluene	72.4% CH ₃ OH, 27.6% C ₆ H ₅ CH ₃	63.7
Methanol–benzene	39.5% CH ₃ OH, 60.5% C ₆ H ₆	58.3
Cyclohexane–ethanol	69.5% C ₆ H ₁₂ , 30.5% C ₂ H ₅ OH	64.9
2-Propanol–water	87.8% (CH ₃) ₂ CHOH, 12.2% H ₂ O	80.4
Butyl acetate–water	72.9% CH ₃ COOC ₄ H ₉ , 27.1% H ₂ O	90.7
Phenol–water	9.2% C ₆ H ₅ OH, 90.8% H ₂ O	99.5

¹ Keep in mind that this distillate is not pure ethanol but is an ethanol–water mixture.

² More examples of azeotropes, with their compositions and boiling points, can be found in the *CRC Handbook of Chemistry and Physics*; also in L. H. Horsley, ed., *Advances in Chemistry Series*, no. 116. Azeotropic Data, III (Washington, DC: American Chemical Society, 1973).

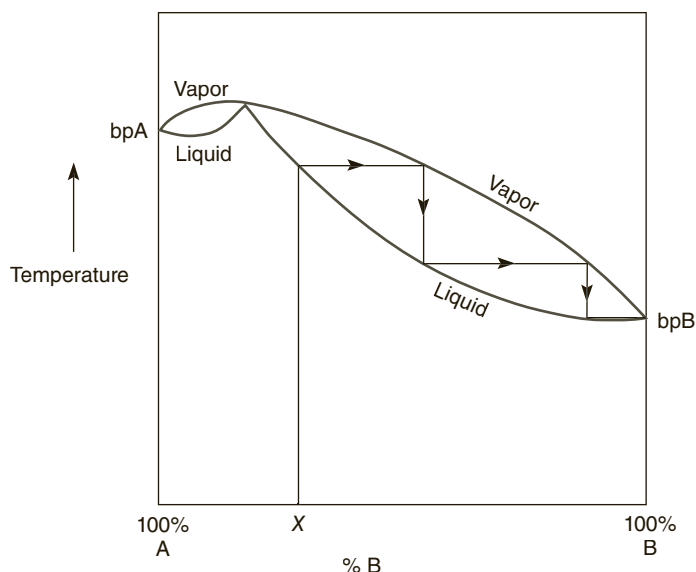


Figure 15.13
A maximum-boiling-point phase diagram.

Table 15.3 Maximum-boiling-point azeotropes

<i>Azeotrope</i>	<i>Composition (Weight Percentage)</i>	<i>Boiling Point (°C)</i>
Acetone–chloroform	20.0% CH ₃ COCH ₃ , 80.0% CHCl ₃	64.7
Chloroform–methyl ethyl ketone	17.0% CHCl ₃ , 83.0% CH ₃ COCH ₂ CH ₃	79.9
Hydrochloric acid	20.2% HCl, 79.8% H ₂ O	108.6
Acetic acid–dioxane	77.0% CH ₃ COCH ₃ , 23.0% C ₄ H ₈ O ₂	119.5
Benzaldehyde–phenol	49.0% C ₆ H ₅ CHO, 51.0% C ₆ H ₅ OH	185.6

expected in the solution. The lower combined vapor pressures cause a higher boiling point than would be characteristic for the components. A two-component maximum-boiling-point azeotrope is illustrated in Figure 15.13. Because the azeotrope has a higher boiling point than any of the components, it will be concentrated in the distillation flask as the distillate (pure B) is removed. The distillation of a solution of composition X would follow to the right along the lines in Figure 15.13. Once the composition of the material remaining in the flask has reached that of the azeotrope, the temperature will rise, and the azeotrope will begin to distill. The azeotrope will continue to distill until all the material in the distillation flask has been exhausted.

Some maximum-boiling-point azeotropes are listed in Table 15.3. They are not nearly as common as minimum-boiling-point azeotropes.

C. Generalizations

There are some generalizations that can be made about azeotropic behavior. They are presented here without explanation, but you should be able to verify them by thinking through each case using the phase diagrams given. (Note that pure A is always to the left of the azeotrope in these diagrams, whereas pure B is to the right of the azeotrope.)

Minimum-boiling-point azeotropes


<i>Initial Composition</i>	<i>Experimental Result</i>
To left of azeotrope	Azeotrope distills first, pure A second
Azeotrope	Not separable
To right of azeotrope	Azeotrope distills first, pure B second

Maximum-boiling-point azeotropes

<i>Initial Composition</i>	<i>Experimental Result</i>
To left of azeotrope	Pure A distills first, azeotrope second
Azeotrope	Not separable
To right of azeotrope	Pure B distills first, azeotrope second

15.8 Azeotropic Distillation: Applications

There are many examples of chemical reactions in which the amount of product is low because of an unfavorable equilibrium. An example is the direct acid-catalyzed esterification of a carboxylic acid with an alcohol:

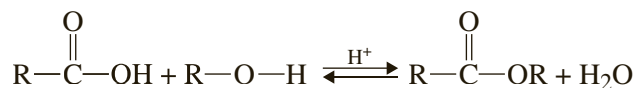

Because the equilibrium does not favor formation of the ester, it must be shifted to the right, in favor of the product, by using an excess of one of the starting materials. In most cases, the alcohol is the least expensive reagent and is the material used in excess. Isopentyl acetate (Experiment 14) is an example of esters prepared by using one of the starting materials in excess.

Another way of shifting the equilibrium to the right is to remove one of the products from the reaction mixture as it is formed. In the preceding example, water can be removed as it is formed by **azeotropic distillation**. A common large-scale method is to use the Dean–Stark water separator shown in Figure 15.14A. In this technique, an inert solvent, commonly benzene or toluene, is added to the reaction mixture contained in the round-bottom flask. The side arm of the water separator is also filled with this solvent. If benzene is used, as the mixture is heated under reflux, the benzene–water azeotrope (bp 69.4°C, Table 15.3) distills out of the flask.³ When the vapor condenses, it enters the side arm directly below the condenser, and water separates from the benzene–water condensate; benzene and water mix as vapors, but they are not miscible as cooled liquids. Once the water (lower phase) separates from the benzene (upper phase), liquid benzene overflows from the side arm back into the flask. The cycle is repeated continuously until no more water forms in the side arm. You may calculate the weight of water that should theoretically be produced and compare this value with the amount of water collected in the side arm. Because the density of water is 1.0, the volume of water collected can be compared directly with the calculated amount, assuming 100% yield.

³ Actually, with ethanol, a lower-boiling-point three-component azeotrope distills at 64.9°C (see Table 15.3). It consists of benzene–water–ethanol. Because some ethanol is lost in the azeotropic distillation, a large excess of ethanol is used in esterification reactions. The excess also helps shift the equilibrium to the right.

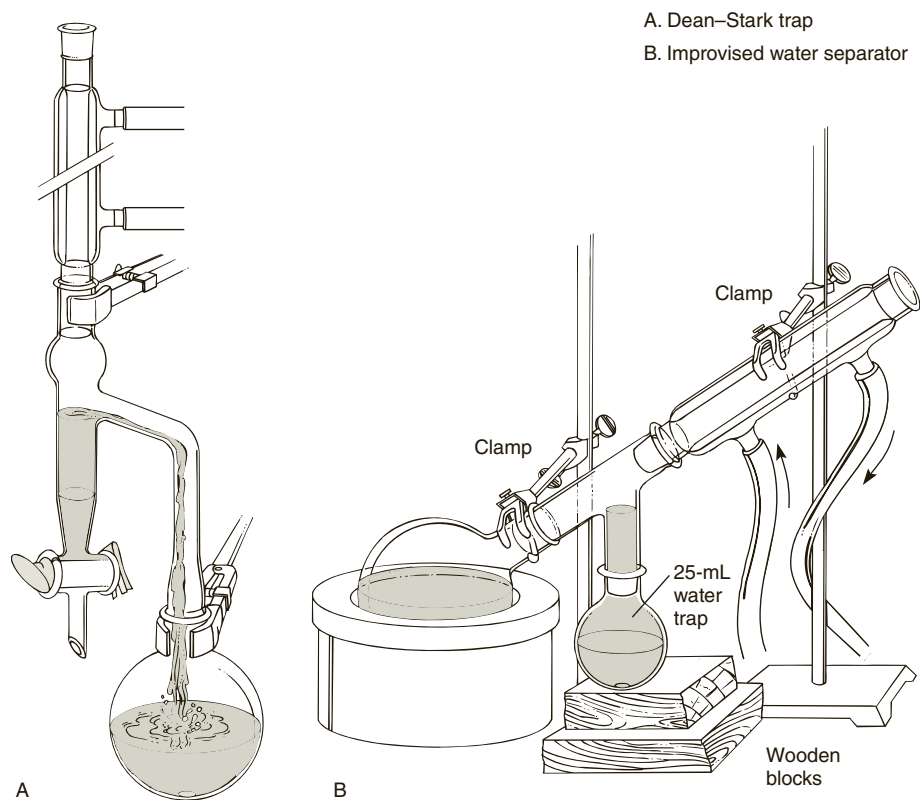


Figure 15.14
Large-scale water separators.

An improvised water separator, constructed from the components found in the traditional organic kit, is shown in Figure 15.14B. Although this requires the condenser to be placed in a nonvertical position, it works quite well.

At the microscale level, water separation can be achieved using a standard distillation assembly with a water condenser and a Hickman head (Figure 15.15). The side-ported variation of the Hickman head is the most convenient one to use for this purpose, but it is not essential. In this variation, you simply remove all the distillate (both solvent and water) several times during the course of the reaction. Use a Pasteur pipette to remove the distillate, as shown in Technique 14, Figure 14.6 Because both the solvent and water are removed in this procedure, it may be desirable to add more solvent from time to time, adding it through the condenser with a Pasteur pipette.

The most important consideration in using azeotropic distillation to prepare an ester (described above) is that the azeotrope containing water must have a **lower boiling point** than the alcohol used. With ethanol, the benzene–water azeotrope boils at a much lower temperature (69.4°C) than ethanol (78.3°C), and the technique previously described works well. With higher-boiling-point alcohols, azeotropic distillation works well because of the large boiling-point difference between the azeotrope and the alcohol.

With methanol (bp 65°C), however, the boiling point of the benzene–water azeotrope is actually *higher* by about 5°C, and methanol distills first. Thus, in esterifications involving methanol, a totally different approach must be taken. For example, you can mix the carboxylic acid, methanol, the acid catalyst, and 1,2-dichloroethane in a conventional reflux apparatus (Technique 7, Figure 7.6) without a water separator. During the reaction, water separates from the


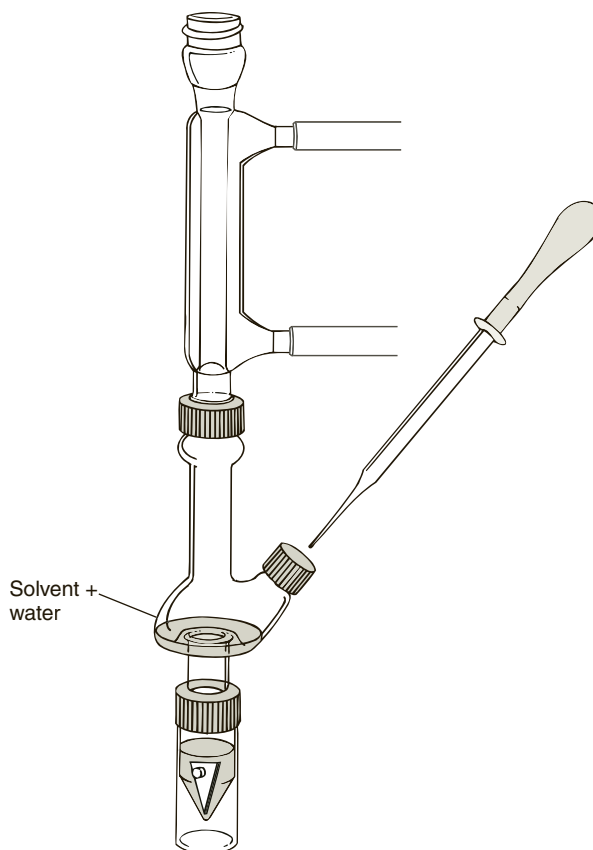


Figure 15.15
Microscale water separator (both layers are removed).

1,2-dichloroethane because it is not miscible; however, the remainder of the components are soluble, so the reaction can continue. The equilibrium is shifted to the right by the “removal” of water from the reaction mixture.

Azeotropic distillation is also used in other types of reactions, such as ketal or acetal formation, and in enamine formation. The use of azeotropic distillation is illustrated in the formation of 2-acetylcyclohexanone via the enamine intermediate. Toluene is used in the azeotropic distillation of water. The Hickman head is used as a water separator.

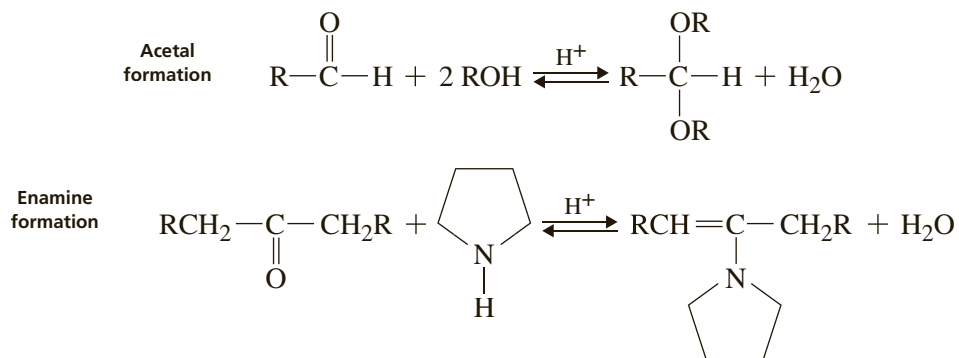

PROBLEMS

1. In the accompanying chart are approximate vapor pressures for benzene and toluene at various temperatures.

<i>Temp</i> (°C)	<i>mm Hg</i>	<i>Temp</i> (°C)	<i>mm Hg</i>
Benzene 30	120	Toluene 30	37
40	180	40	60
50	270	50	95
60	390	60	140
70	550	70	200
80	760	80	290
90	1010	90	405
100	1340	100	560
		110	760

- What is the mole fraction of each component if 3.9 g of benzene C_6H_6 is dissolved in 4.6 g of toluene C_7H_8 ?
 - Assuming that this mixture is ideal—that is, it follows Raoult's Law—what is the partial vapor pressure of benzene in this mixture at 50°C?
 - Estimate to the nearest degree the temperature at which the vapor pressure of the solution equals 1 atm (bp of the solution).
 - Calculate the composition of the vapor (mole fraction of each component) that is in equilibrium in the solution at the boiling point of this solution.
 - Calculate the composition in weight percentage of the vapor that is in equilibrium with the solution.
- Estimate how many theoretical plates are needed to separate a mixture that has a mole fraction of B equal to 0.70 (70% B) in Figure 15.4.
 - Two moles of sucrose are dissolved in 8 moles of water. Assume that the solution follows Raoult's Law and that the vapor pressure of sucrose is negligible. The boiling point of water is 100°C. The distillation is carried out at 1 atm (760 mm Hg).
 - Calculate the vapor pressure of the solution when the temperature reaches 100°C.
 - What temperature would be observed during the entire distillation?
 - What would be the composition of the distillate?
 - If a thermometer were immersed below the surface of the liquid of the boiling flask, what temperature would be observed?
 - Explain why the boiling point of a two-component mixture rises slowly throughout a simple distillation when the boiling-point differences are not large.
 - Given the boiling points of several known mixtures of A and B (mole fractions are known) and the vapor pressures of A and B in the pure state (P_A° and P_B°) at these same temperatures, how would you construct a boiling-point-composition phase diagram for A and B? Give a stepwise explanation.

6. Describe the behavior on distillation of a 98% ethanol solution through an efficient column. Refer to Figure 15.12.
7. Construct an approximate boiling-point-composition diagram for a benzene-methanol system. The mixture shows azeotropic behavior (see Table 15.3). Include on the graph the boiling points of pure benzene and pure methanol and the boiling point of the azeotrope. Describe the behavior on distillation of a mixture that is initially rich in benzene (90%) and then for a mixture that is initially rich in methanol (90%).
8. Construct an approximate boiling-point-composition diagram for an acetone-chloroform system, which forms a maximum boiling azeotrope (Table 15.4). Describe the behavior on distillation of a mixture that is initially rich in acetone (90%), and then describe the behavior of a mixture that is initially rich in chloroform (90%).
9. Two compounds have boiling points of 130 and 150°C. Estimate the number of theoretical plates needed to separate these substances in a fractional distillation.
10. A spinning-band column has an HETP of 0.25 in./plate. If the column has 12 theoretical plates, how long is it?

Vacuum Distillation, Manometers

Vacuum distillation (distillation at reduced pressure) is used for compounds that have high boiling points (above 200°C). Such compounds often undergo thermal decomposition at the temperatures required for their distillation at atmospheric pressure. The boiling point of a compound is lowered substantially by reducing the applied pressure. Vacuum distillation is also used for compounds that, when heated, might react with the oxygen present in air. It is also used when it is more convenient to distill at a lower temperature because of experimental limitations. For instance, a heating device may have difficulty heating to a temperature in excess of 250°C.

The effect of pressure on the boiling point is discussed more thoroughly in Technique 13, Section 13.1. A nomograph is given (Figure 13.2) that allows you to estimate the boiling point of a liquid at a pressure different from the one at which it is reported. For example, a liquid reported to boil at 200°C at 760 mm Hg would be expected to boil at 90°C at 20 mm Hg. This is a significant decrease in temperature, and it would be advantageous to use a vacuum distillation if any problems were to be expected. Counterbalancing this advantage, however, is the fact that separations of liquids of different boiling points may not be as effective with a vacuum distillation as with a simple distillation.

16.1 Microscale Methods

When working with glassware that is to be evacuated, you should wear safety glasses at all times. There is always danger of an implosion.

CAUTION



Safety glasses must be worn at all times during vacuum distillation.

It is a good idea to work in a hood when performing a vacuum distillation. If the experiment will involve high temperatures (> 220°C) for distillation or an extremely low pressure (< 0.1 mm Hg), for your own safety you should definitely work in a hood, behind a shield.

A basic apparatus similar to the one shown in Figure 16.1 may be used for microscale vacuum distillations. As is the case for simple distillation, this apparatus uses the Hickman head as a means to reduce the length of the vapor path. The major difference to be found when comparing this assembly to one for simple distillation (Technique 14, Figure 14.8) is that the opening to the atmosphere has been replaced by a connection to a vacuum source (top right-hand side). The usual sources of vacuum are the aspirator (Technique 8, Section 8.5), a mechanical vacuum pump, or a “house” vacuum line (one piped directly to the laboratory bench). The aspirator is probably the simplest of these sources and the vacuum source most likely to be available. However, if pressures below 10–20 mm Hg are required, a vacuum pump must be used.

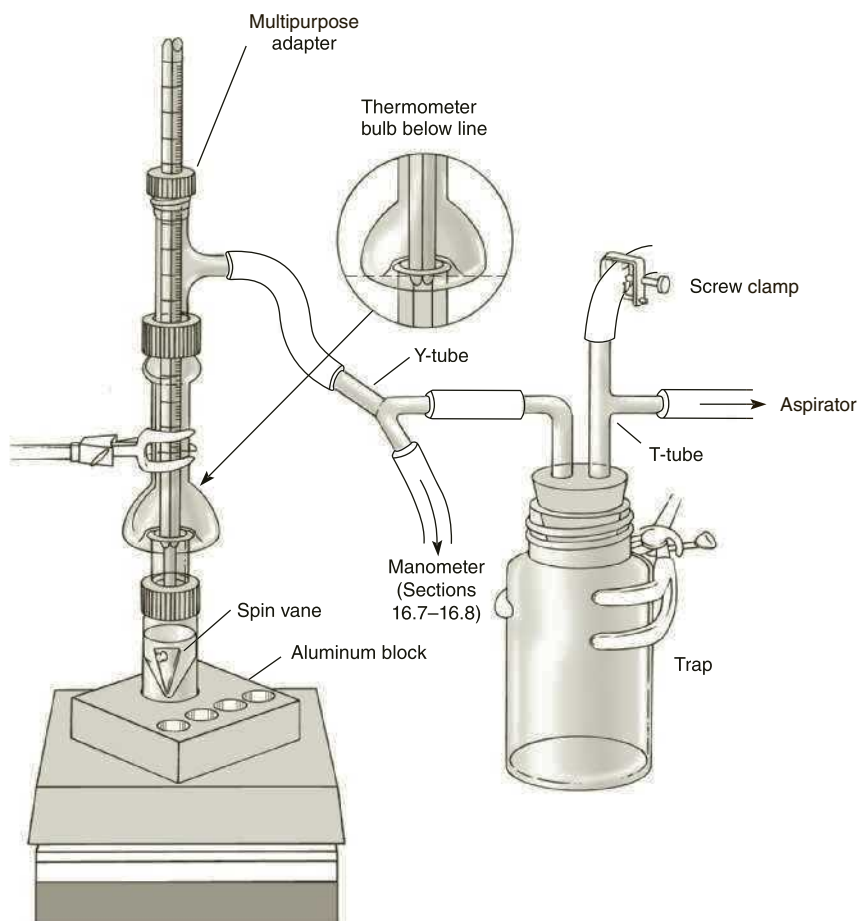


Figure 16.1
Reduced-pressure microscale distillation (internal monitoring of temperature).

Assembling the Apparatus. When assembling an apparatus for vacuum distillation, it is important that all joints and connections be airtight. The joints in the newest microscale kits are standard-taper ground-glass joints, with a compression cap that contains an O-ring seal. Glassware that contains this type of compression joint will hold a vacuum quite easily. Under normal conditions, it is not necessary to grease these joints.

NOTE: Normally, you should not grease joints. It is necessary to grease the joints in a vacuum distillation only if you cannot achieve the desired pressure without using grease.

If you must grease joints, take care not to use too much grease. You are working with small quantities of liquid in a microscale distillation, and the grease can become a very serious contaminant if it oozes out the bottom of the joints into your system. Apply a small amount of grease (thin film) completely around the top of the *inner* joint; then mate the joints and turn them slightly to spread the grease evenly. If you have used the correct amount of grease, it will not ooze out the bottom; rather, the entire joint will appear clear and without striations or uncovered areas.

Make doubly sure that any connections to pressure tubing are tight. The pressure tubing itself should be relatively new and without cracks. If the tubing shows cracks when you stretch or bend it, it may be old and leak air into the system. Glass

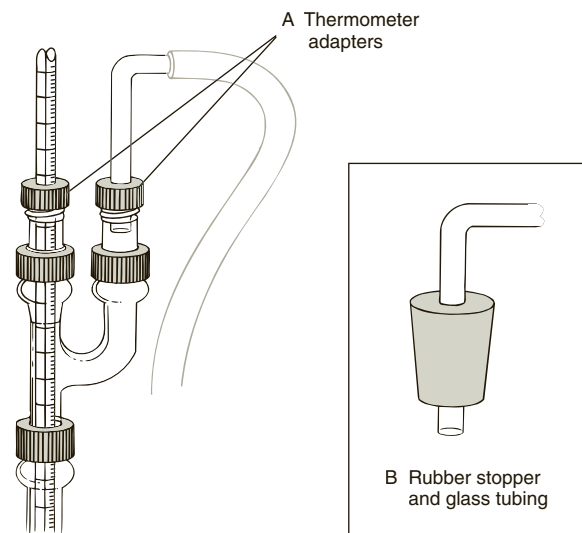


Figure 16.2
Alternative vacuum connections.

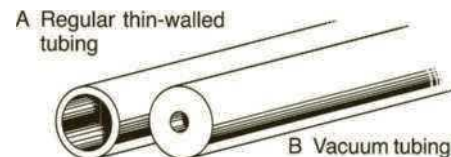


Figure 16.3
Comparison of tubing.

tubing should fit securely into any rubber stoppers. If you can move the tubing up and down with only gentle force, it is too loose, and you should obtain a larger size. Check all glassware to be sure there are no cracks and that there are no chips in the standard-taper joints. Cracked glassware may break when evacuated.

Connecting to Vacuum. In Figure 16.1, the connection to vacuum has been made using a multipurpose adapter (see Technique 14, Figure 14.9B). If a multipurpose adapter is not available, an alternative method uses a Claisen head and two thermometer adapters (Figure 16.2). If two thermometer adapters are not available, a #0 rubber stopper fitted with glass tubing can be used.

Whichever is used, the connection to the vacuum source is made using **pressure tubing**. Pressure tubing (also called vacuum tubing), unlike the more common thin-walled tubing used to carry water or gas, has heavy walls that will not collapse inward when it is evacuated. Compare the two types of tubing shown in Figure 16.3.

Water Trap. If an aspirator is used as a source of vacuum, a water trap must be placed between it and the distillation assembly. A commonly used type of water trap is shown at the bottom right of Figure 16.1. Variations in water pressure are to be expected when using an aspirator. If the pressure drops low enough, the vacuum in the system will draw water from the aspirator into the connecting line. The trap allows you to see this happening and take corrective action (prevent water from entering the distillation apparatus). The correct action for anything but a small amount of water is to “vent the system.” This can be accomplished by opening the screw clamp at the top of the trap to let air into the system. When performing a vacuum distillation, you should also realize that the system should always be vented before stopping the aspirator. If you turn off the aspirator while the system is still under vacuum, water will be drawn into the connecting line and trap.

Manometer Connection. A Y-tube is shown in the line from the apparatus to the trap. This branching connection is optional but is required if you wish to monitor

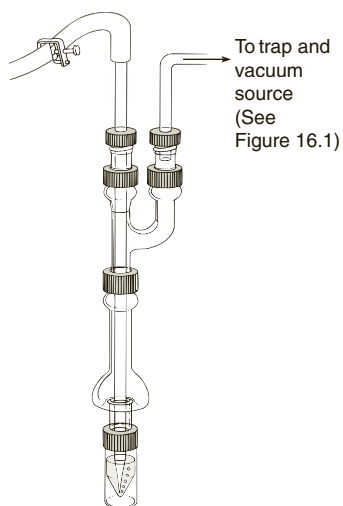


Figure 16.4
Use of ebulliator tube
instead of thermometer.

the actual pressure of the system using a manometer. The operation of manometers is discussed in Sections 16.7 and 16.8.

Thermometer Placement. If a thermometer is used, be sure that the bulb is placed in the stem of the Hickman head just below the well. If it is placed higher, it may not be surrounded by a constant stream of vapor from the material being distilled. If the thermometer is not exposed to a continuous stream of vapor, it may not reach temperature equilibrium. As a result, the temperature reading would be incorrect (low).

Preventing Bumpover. When a distillation flask is heated, there is always the possibility that the boiling action will become too vigorous (mainly due to superheating) and “bump” some of the undistilled liquid up into the Hickman head. The simplest way to prevent bumping is to stir the boiling liquid with a magnetic spin vane. Stirring rapidly will distribute the heat evenly, keep the boiling action smooth, and prevent bumping. Boiling stones cannot be used for this purpose in a vacuum distillation; they do not work in vacuum. In a conventional vacuum distillation (macroscale), it is customary to maintain smooth boiling action by using an **ebulliator tube**.

The ebulliator tube agitates the boiling solution by providing a small, continuous stream of air bubbles. Figure 16.4 shows how a microscale vacuum distillation may be modified to use an ebulliator tube. The amount of air (rate of bubbles) provided by the ebulliator is adjusted by either tightening or loosening the screw clamp at the top. A Pasteur pipette makes an excellent ebulliator tube. As Figure 16.4 shows, the ebulliator tube replaces the thermometer. Hence, the ebulliator should be used only when internal monitoring of temperature is not required. In practice, although this method works satisfactorily, better results are obtained by stirring and distilling slowly.

NOTE: Heating slowly helps to avoid bumping.

16.2 Semi-microscale and Macroscale Equipment

A vacuum distillation apparatus using the components of the traditional organic laboratory kit is shown in Figure 16.5. It uses the ebulliator tube, the Claisen head, and a thermometer for internal temperature monitoring. A water condenser is shown, but with high-boiling liquids, this apparatus may be simplified by removing the water condenser. A special vacuum adapter allows connection to the manometer and vacuum source.

The Claisen head is used in larger-scale vacuum distillations because it allows the use of an ebulliator tube. The bend it provides in the distilling path helps to prevent bumpover. Because the Claisen head increases the holdup of the system, it cannot be used with very-small-scale distillations (<10 mL).

16.3 Stepwise Instructions for Microscale Vacuum Distillations

The following set of instructions is a step-by-step account of how to carry out a vacuum distillation. The microscale apparatus illustrated in Figure 16.1 will be used; however, the procedures apply to any vacuum distillation.

CAUTION

Safety glasses must be worn at all times during vacuum distillation.



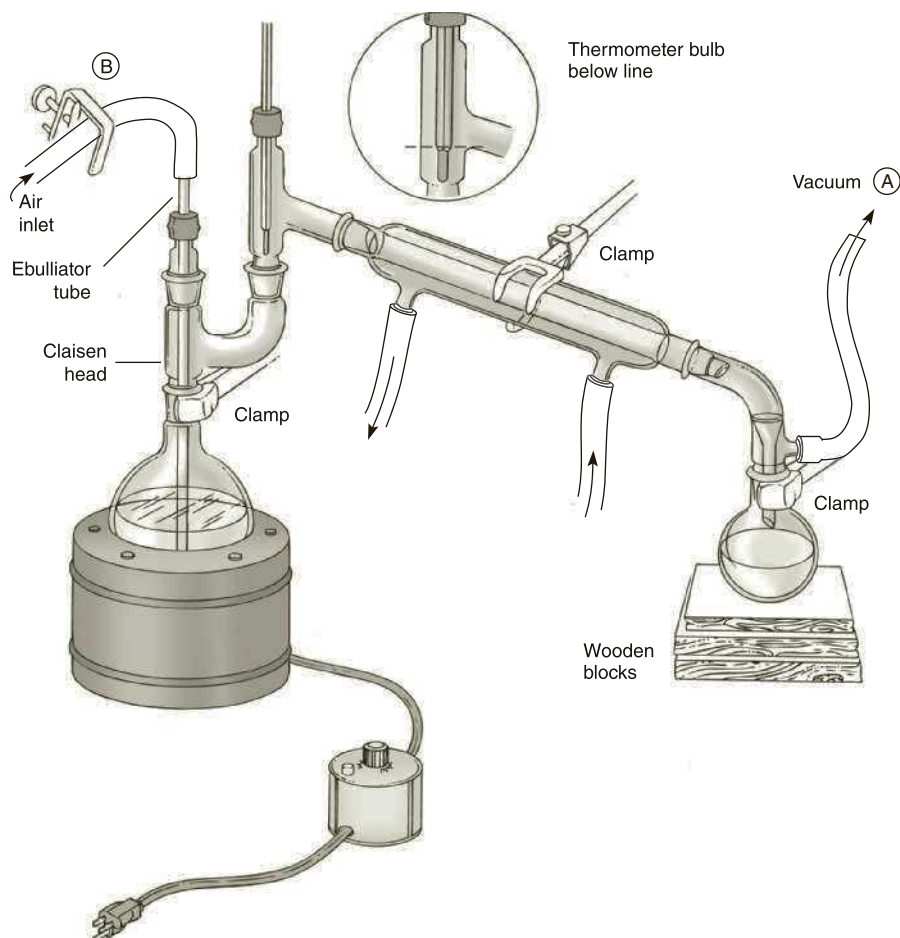


Figure 16.5
 Macroscale vacuum distillation using the standard organic laboratory kit.

A. Evacuating the Apparatus

1. Assemble the apparatus as shown in Figure 16.1. It should be held with a clamp attached to the top of the Hickman head and placed *above* the aluminum block.

NOTE: If you expect the temperature of the distillation to rise above 150°C, omit the threaded cap and O-ring between the conical vial and Hickman head. They will melt at high temperature.

2. If the sample contains solvent, concentrate the sample to be distilled in the conical vial (or round-bottom flask) that you are using. Use one of the solvent-removal methods discussed in Technique 7. If you have a large volume of solvent to evaporate and the sample does not fit in the conical vial, you must use an Erlenmeyer flask first and then transfer the sample to the conical vial. (Be sure to rinse the Erlenmeyer flask with a little solvent and then reevaporate in the conical vial.) As a rule, the distillation vial or flask should be no more than two-thirds full.
3. Attach the conical vial (or flask) to the apparatus and make sure all joints are sealed.
4. Turn the aspirator on to the maximum extent.

5. Close the screw clamp on the water trap very tightly. (If you are using an ebulliator tube as in Figure 16.4, you next regulate the rate of bubbling by adjusting the tightness of the screw clamp at the top of the tube.)
6. Using the manometer, observe the pressure. It may take a few minutes to remove any residual solvent and evacuate the system. If the pressure is not satisfactory, check all connections to see whether they are tight. (Readjust the ebulliator tube if necessary.)

NOTE: Do not proceed until you have a good vacuum.

B. Beginning Distillation

7. Lower the apparatus into the aluminum block and begin to heat. Place the external thermometer in the block now if you wish.

CAUTION



You should not use a *mercury* thermometer with an aluminum block. If it breaks, the mercury will vaporize on the hot surface. Instead, use a nonmercury glass thermometer, a metal dial thermometer, or a digital electronic temperature-measuring device.

8. Increase the temperature of the heat source until you begin to see distillate collect in the well of the Hickman head. (Observe very carefully; liquid may appear almost “magically” without any sign of boiling or any obvious reflux ring.)
9. If you are using a thermometer, record the temperature and pressure when distillate begins to appear. (If you are not using an internal thermometer, record the external temperature. If you have two thermometers, record both temperatures.)

C. Collecting a Fraction

10. To collect a fraction, raise the apparatus above the aluminum block and allow it to cool a bit before opening it.
11. Open the screw clamp on the water trap to allow air to enter the system. (If you are using an ebulliator tube, you also need to open the screw clamp at its *top immediately*, or the liquid in the distillation flask will be forced upward into it.)
12. Partially disassemble the apparatus and remove the fraction with a Pasteur pipette, as shown in Technique 14, Figure 14.6A. (If you have a Hickman head with a side port, you may simply open the side port to remove the fraction. This is shown in Figure 14.6C.)

NOTE: If you do not intend to collect a second fraction, go directly to steps 18–20.

13. Reassemble the apparatus (or close the side port) and tighten the clamp at the top of the ebulliator tube.
14. Tighten the screw clamp on the water trap and reestablish the desired pressure. If the pressure is not satisfactory, check all connections to make sure they are sealed.

- Lower the apparatus back into the aluminum block and continue the distillation.

D. Shutdown

- At the end of the distillation, raise the apparatus from the aluminum block and allow it to cool. Also let the aluminum block cool.
- Open the screw clamp on the water trap first, and then immediately open the one at the top of the ebulliator tube.
- Turn off the water at the aspirator. (Do not do this before step 17!)
- Remove any distilled material by one of the methods shown in Figure 14.6.
- Disassemble the apparatus and clean all glassware as soon as possible to prevent the joints from sticking.

NOTE: If you used grease, thoroughly clean all grease off the joints, or it will contaminate your samples in other procedures.

16.4 Rotary Fraction Collectors

With the types of apparatus we have discussed previously, the vacuum must be stopped to remove fractions when a new substance (fraction) begins to distill. Quite a few steps are required to perform this change, and it is quite inconvenient when there are several fractions to be collected. Two pieces of semimicroscale apparatus that are designed to alleviate the difficulty of collecting fractions while working under vacuum are shown in Figure 16.6. The collector, which is shown to the right, is sometimes called a “cow” because of its appearance. With these rotary fraction collecting devices, all you need to do is rotate the device to collect fractions.

16.5 Bulb-to-Bulb Distillation

The ultimate in microscale methods is to use a bulb-to-bulb distillation apparatus. This apparatus is shown in Figure 16.7. The sample to be distilled is placed in the

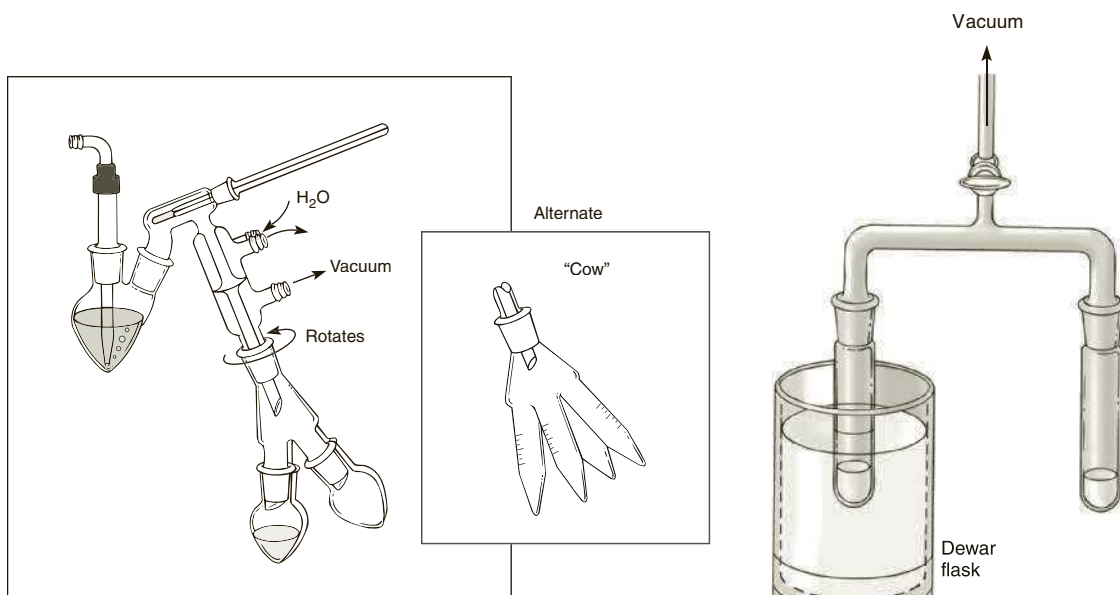


Figure 16.6
Rotary fraction collector.

Figure 16.7
Bulb-to-bulb distillation.

glass container attached to one of the arms of the apparatus. The sample is frozen solid, usually by using liquid nitrogen, but dry ice in 2-propanol or an ice–salt-water mixture may also be used. The coolant container shown in the figure is a **Dewar flask**. The Dewar flask has a double wall with the space between the walls evacuated and sealed. A vacuum is a very good thermal insulator, and there is little heat loss from the cooling solution.

After the sample is frozen, the entire apparatus is evacuated by opening the stopcock. When the evacuation is complete, the stopcock is closed, and the Dewar flask is removed. The sample is allowed to thaw, and then it is frozen again. This freeze–thaw–freeze cycle removes any air or gases that were trapped in the frozen sample. Next, the stopcock is opened to evacuate the system again. When the second evacuation is complete, the stopcock is closed, and the Dewar flask is moved to the other arm to cool the empty container. As the sample warms, it will vaporize, travel to the other side, and be frozen or liquefied by the cooling solution. This transfer of the liquid from one arm to the other may take quite a while, but no heating is required.

The bulb-to-bulb distillation is most effective when liquid nitrogen is used as coolant and when the vacuum system can achieve a pressure of 10^{-3} mm Hg or lower. This requires a vacuum pump; an aspirator cannot be used.

16.6 The Mechanical Vacuum Pump

The aspirator is not capable of yielding pressures below about 5 mm Hg. This is the vapor pressure of water at 0°C, and water freezes at this temperature. A more realistic value of pressure for an aspirator is about 20 mm Hg. When pressures below 20 mm Hg are required, a vacuum pump will have to be employed. Figure 16.8 illustrates a mechanical vacuum pump and its associated glassware. The vacuum pump operates on a principle similar to that of the aspirator, but the vacuum pump uses a high-boiling oil, rather than water, to remove air from the attached system. The oil used in a vacuum pump, a silicone oil or a high-molecular-weight hydrocarbon-based oil, has

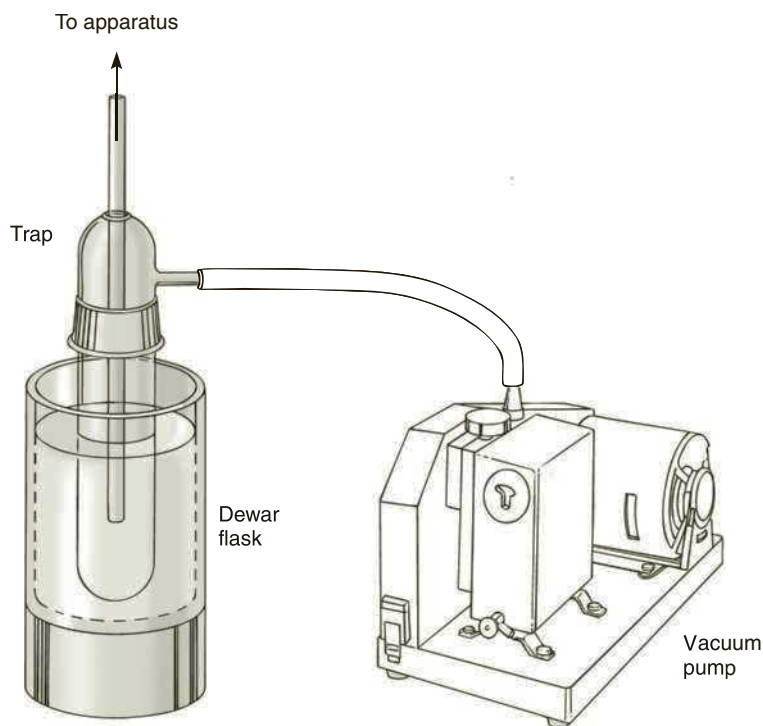


Figure 16.8
A vacuum pump and its trap.

a very low vapor pressure, and very low system pressures can be achieved. A good vacuum pump, with new oil, can achieve pressures of 10^{-3} or 10^{-4} mm Hg. Instead of the oil being discarded as it is used, it is recycled continuously through the system.

A cooled trap is required when using a vacuum pump. This trap protects the oil in the pump from any vapors that may be present in the system. If vapors from organic solvents, or from the organic compounds being distilled, dissolve in the oil, the vapor pressure of the oil will increase, rendering it less effective. A special type of vacuum trap is illustrated in Figure 16.8. It is designed to fit into an insulated Dewar flask so that the coolant will last for a long period. At a minimum, this flask should be filled with ice water, but a dry ice–acetone mixture or liquid nitrogen is required to achieve lower temperatures and better protect the oil. Often, two traps are used; the first trap contains ice water and the second trap dry ice–acetone or liquid nitrogen. The first trap liquefies low-boiling vapors that might freeze or solidify in the second trap and block it.

16.7 The Closed-End Manometer

The principal device used to measure pressures in a vacuum distillation is the **closed-end manometer**. Two basic types are shown in Figures 16.9 and 16.10. The manometer shown in Figure 16.9 is widely used because it is relatively easy to construct. It consists of a U-tube that is closed at one end and mounted on a wooden support. You can construct the manometer from 9-mm glass capillary tubing and fill it, as shown in Figure 16.11.

CAUTION



Mercury is a very toxic metal with cumulative effects. Because mercury has a high vapor pressure, it must not be spilled in the laboratory. You must not touch it. Seek immediate help from an instructor in case of a spill or if you break a manometer. Spills must be cleaned immediately.

A small filling device is connected to the U-tube with pressure tubing. The U-tube is evacuated with a good vacuum pump; then the mercury is introduced by tilting the mercury reservoir.

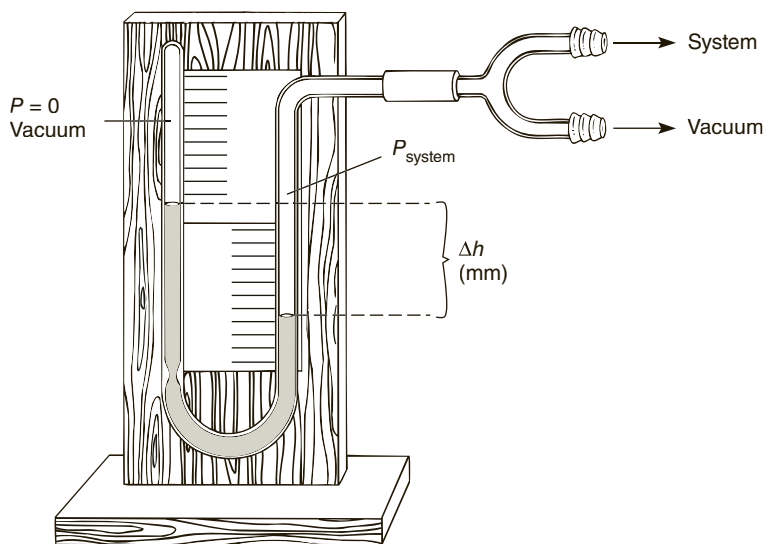


Figure 16.9
A simple U-tube manometer.

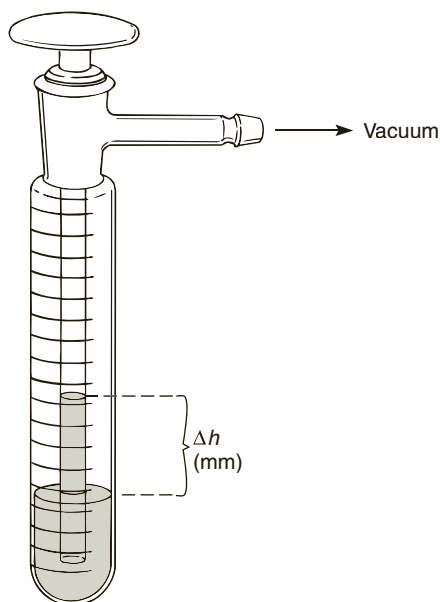


Figure 16.10
Commercial “stick” manometer.

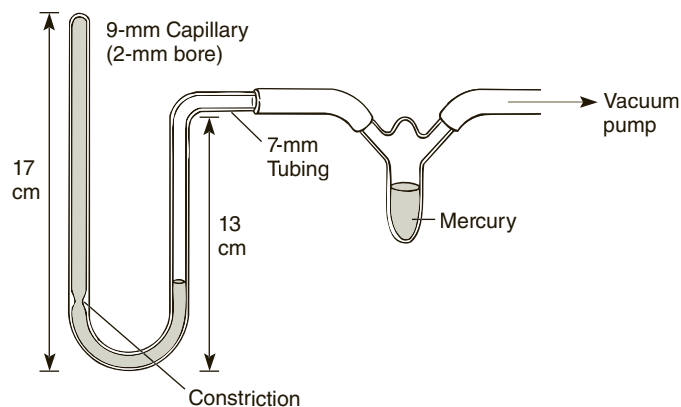


Figure 16.11
Filling a U-tube manometer.

NOTE: The entire filling operation should be conducted in a shallow pan in order to contain any spills that might occur.

Enough mercury should be added to form a column about 20 cm in total length. When the vacuum is interrupted by admitting air, the mercury is forced by atmospheric pressure to the end of the evacuated tube. The manometer is then ready for use. The constriction shown in Figure 16.11 helps to protect the manometer against breakage when the pressure is released. Be sure that the column of mercury is long enough to pass through this constriction.

When an aspirator or any other vacuum source is used, a manometer can be connected into the system. As the pressure is lowered, the mercury rises in the right tube and drops in the left tube until Δh corresponds to the approximate pressure of the system (see Figure 16.9).

$$\Delta h = (P_{\text{system}} - P_{\text{reference arm}}) = (P_{\text{system}} - 10^{-3} \text{ mm Hg}) \approx P_{\text{system}}$$

A short piece of metric ruler or a piece of graph paper ruled in millimeter squares is mounted on the support board to allow Δh to be read. No addition or subtraction is necessary, because the reference pressure (created by the initial evacuation when filling) is approximately zero (10^{-3} mm Hg) when referred to readings in the 10- to 50-mm Hg range. To determine the pressure, count the number of millimeter squares beginning at the top of the mercury column on the left and continuing downward to the top of the mercury column on the right. This is the height difference Δh , and it gives the pressure in the system directly.

A commercial counterpart to the U-tube manometer is shown in Figure 16.10. With this manometer, the pressure is given by the difference in the mercury levels in the inner and outer tubes.

The manometers described here have a range of about 1–150 mm Hg in pressure. They are convenient to use when an aspirator is the source of vacuum. For high-vacuum systems (pressures below 1 mm Hg), a more elaborate manometer or an electronic measuring device must be used. These devices will not be discussed here.

16.8 Connecting and Using a Manometer

The most common use of a closed-end manometer is to monitor pressure during a reduced-pressure distillation. The manometer is placed in a vacuum distillation system, as shown in Figure 16.12. Generally, an aspirator is the source of vacuum. Both the manometer and the distillation apparatus should be protected by a trap from possible backups in the water line. Alternatives to the trap arrangements shown in Figure 16.12 appear in Figure 16.1. Notice in each case that the trap has a device (screw clamp or stopcock) for opening the system to the atmosphere. This is especially important in using a manometer because you should always make pressure changes slowly. If this is not done, there is a danger of spraying mercury throughout the system, breaking the manometer, or spurring mercury into the room. In the closed-end manometer, if the system is opened suddenly, the mercury rushes to the closed end of the U-tube. The mercury rushes with such speed and force that the end will be broken out of the manometer. Air should be admitted *slowly* by opening the valve cautiously. In a similar fashion, the valve should be closed slowly when the vacuum is being started, or mercury may be forcefully drawn into the system through the open end of the manometer.

If the pressure in a reduced-pressure distillation is lower than that desired, it is possible to adjust it by means of a **bleed valve**. The stopcock can serve this function in Figure 16.12 if it is opened only a small amount. In those systems with a screw clamp on the trap (Figure 16.1), remove the screw clamp from the trap valve and

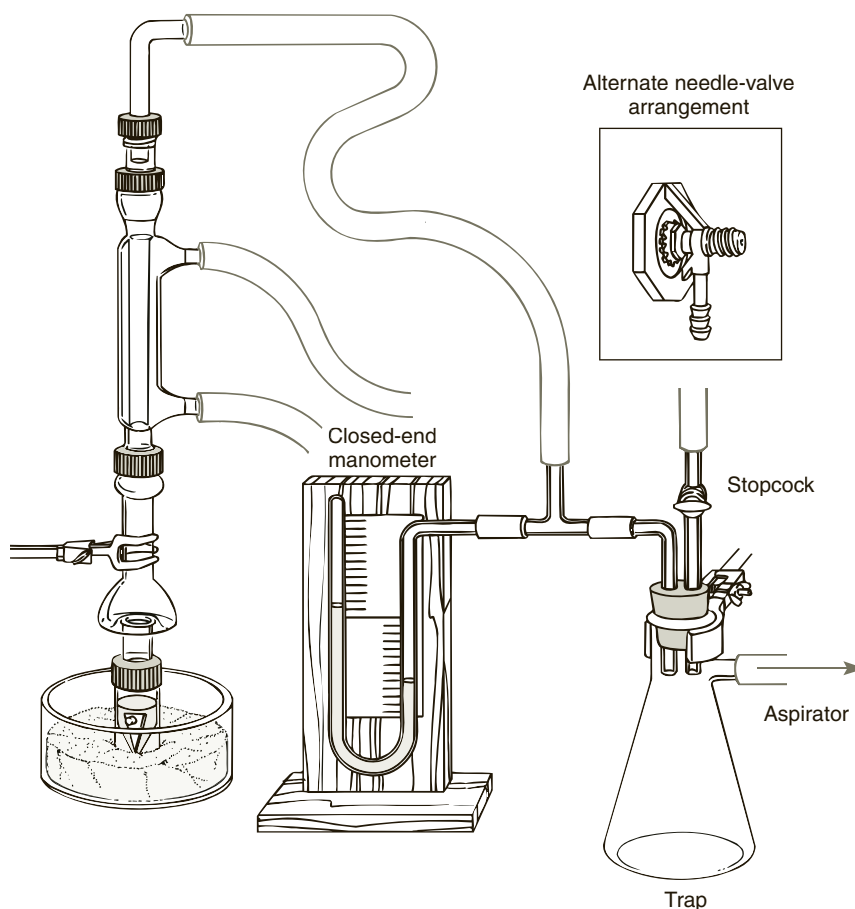


Figure 16.12

Connecting a manometer to the system. To construct a “bleed,” the needle valve may replace the stopcock.

attach the base of a Tirrill-style Bunsen burner. The needle valve in the base of the burner can be used to adjust precisely the amount of air that is admitted (bled) to the system and hence control the pressure.

PROBLEMS

1. Give some reasons that would lead you to purify a liquid by using vacuum distillation rather than by using simple distillation.
2. When using an aspirator as a source of vacuum in a vacuum distillation, do you turn off the aspirator before venting the system? Explain.
3. A compound was distilled at atmospheric pressure and had a boiling range of 310–325°C. What would be the approximate boiling range of this liquid if it were distilled under vacuum at 20 mm Hg?
4. Boiling stones generally do not work when a vacuum distillation is performed. What substitutes may be used?
5. What is the purpose of the trap that is used during a vacuum distillation performed with an aspirator?

Sublimation

In Technique 13, the influence of temperature on the change in vapor pressure of a liquid was considered (see Figure 13.1). It was shown that the vapor pressure of a liquid increases with temperature. Because the boiling point of a liquid occurs when its vapor pressure is equal to the applied pressure (normally atmospheric pressure), the vapor pressure of a liquid equals 760 mm Hg at its boiling point. The vapor pressure of a solid also varies with temperature. Because of this behavior, some solids can pass directly into the vapor phase without going through a liquid phase. This process is called **sublimation**. Because the vapor can be resolidified, the overall vaporization–solidification cycle can be used as a purification method. The purification can be successful only if the impurities have significantly lower vapor pressures than the material being sublimed.

17.1 Vapor-Pressure Behavior of Solids and Liquids

In Figure 17.1, vapor-pressure curves for solid and liquid phases for two different substances are shown. Along lines *AB* and *DF*, the sublimation curves, the solid and vapor are at equilibrium. To the left of these lines, the solid phase exists, and

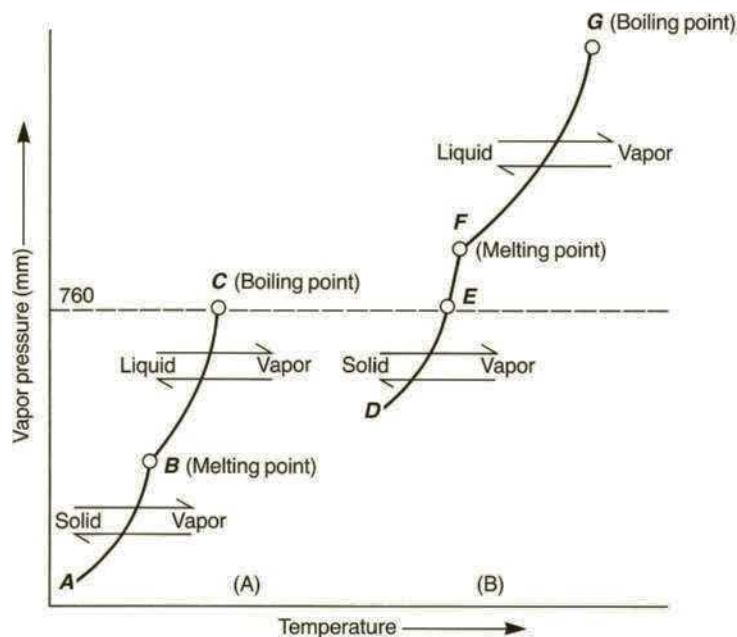


Figure 17.1

Vapor-pressure curves for solids and liquids. (A) Substance shows normal solid to liquid to gas transitions at 760 mm Hg pressure. (B) Substance shows a solid to gas transition at 760 mm Hg pressure.

to the right of these lines, the vapor phase is present. Along lines *BC* and *FG*, the liquid and vapor are at equilibrium. To the left of these lines, the liquid phase exists, and to the right, the vapor is present. The two substances vary greatly in their physical properties, as shown in Figure 17.1.

In the first case (Figure 17.1A), the substance shows normal change-of-state behavior on being heated, going from solid to liquid to gas. The dashed line, which represents an atmospheric pressure of 760 mm Hg, is located *above* the melting point *B* in Figure 17.1A. Thus, the applied pressure (760 mm Hg) is *greater* than the vapor pressure of the solid–liquid phase at the melting point. Starting at *A*, as the temperature of the solid is raised, the vapor pressure increases along *AB* until the solid is observed to melt at *B*. At *B* the vapor pressures of *both* the solid and liquid are identical. As the temperature continues to rise, the vapor pressure will increase along *BC* until the liquid is observed to boil at *C*. The description given is for the “normal” behavior expected for a solid substance. All three states (solid, liquid, and gas) are observed sequentially during the change in temperature.

In the second case (Figure 17.1B), the substance develops enough vapor pressure to vaporize completely at a temperature below its melting point. The substance shows a solid-to-gas transition only. The dashed line is now located *below* the melting point *F* of this substance. Thus, the applied pressure (760 mm Hg) is *less* than the vapor pressure of the solid–liquid phase at the melting point. Starting at *D*, the vapor pressure of the solid rises as the temperature increases along line *DF*. However, the vapor pressure of the solid reaches atmospheric pressure (point *E*) *before* the melting point at *F* is attained. Therefore, sublimation occurs at *E*. No melting behavior will be observed at atmospheric pressure for this substance. For a melting point to be reached and the behavior along line *FG* to be observed, an applied pressure greater than the vapor pressure of the substance at point *F* would be required. This could be achieved by using a sealed pressure apparatus.

The sublimation behavior just described is relatively rare for substances at atmospheric pressure. Several compounds exhibiting this behavior—carbon dioxide, perfluorocyclohexane, and hexachloroethane—are listed in Table 17.1. Notice that these compounds have vapor pressures *above* 760 mm Hg at their melting points. In other words, their vapor pressures reach 760 mm Hg below their melting points, and they sublime rather than melt. Anyone trying to determine the melting point of hexachloroethane at atmospheric pressure will see vapor pouring from the end of the melting-point tube! With a sealed capillary tube, the melting point of 186°C is observed.

Table 17.1 Vapor pressures of solids at their melting points

<i>Compound</i>	<i>Vapor Pressure of Solid at MP (mm Hg)</i>	<i>Melting Point (°C)</i>
Carbon dioxide	3876 (5.1 atm)	−57
Perfluorocyclohexane	950	59
Hexachloroethane	780	186
Camphor	370	179
Iodine	90	114
Naphthalene	7	80
Benzoic acid	6	122
<i>p</i> -Nitrobenzaldehyde	0.009	106

17.2 Sublimation Behavior of Solids

Sublimation is usually a property of relatively nonpolar substances that also have highly symmetrical structures. Symmetrical compounds have relatively high melting points and high vapor pressures. The ease with which a substance can escape from the solid state is determined by the strength of intermolecular forces. Symmetrical molecular structures have a relatively uniform distribution of electron density and a small dipole moment. A smaller dipole moment means a higher vapor pressure because of lower electrostatic attractive forces in the crystal.

Solids sublime if their vapor pressures are greater than atmospheric pressure at their melting points. Some compounds with the vapor pressures at their melting points are listed in Table 17.1. The first three entries in the table were discussed in Section 17.1. At atmospheric pressure they would sublime rather than melt, as shown in Figure 17.1B.

The next four entries in Table 17.1 (camphor, iodine, naphthalene, and benzoic acid) exhibit typical change-of-state behavior (solid, liquid, and gas) at atmospheric pressure, as shown in Figure 17.1A. These compounds sublime readily under reduced pressure, however. Vacuum sublimation is discussed in Section 17.3.

Compared with many other organic compounds, camphor, iodine, and naphthalene have relatively high vapor pressures at relatively low temperatures. For example, they have a vapor pressure of 1 mm Hg at 42, 39, and 53°C, respectively. Although this vapor pressure does not seem very large, it is high enough to lead, after a time, to **evaporation** of the solid from an open container. Mothballs (naphthalene and 1,4-dichlorobenzene) show this behavior. When iodine stands in a closed container over a period of time, you can observe movement of crystals from one part of the container to another.

Although chemists often refer to any solid–vapor transition as sublimation, the process just described for camphor, iodine, and naphthalene is really an **evaporation** of a solid. Strictly speaking, a sublimation point is like a melting point or a boiling point. It is defined as the point at which the vapor pressure of the solid *equals* the applied pressure. Many liquids readily evaporate at temperatures far below their boiling points. It is, however, much less common for solids to evaporate. Solids that readily sublime (evaporate) must be stored in sealed containers. When the melting point of such a solid is being determined, some of the solid may sublime and collect toward the open end of the melting-point tube while the rest of the sample melts. To solve the sublimation problem, one seals the capillary tube or rapidly determines the melting point. It is possible to use the sublimation behavior to purify a substance. For example, at atmospheric pressure, camphor can be readily sublimed, just below its melting point at 175°C. At 175°C the vapor pressure of camphor is 320 mm Hg. The vapor solidifies on a cool surface.

17.3 Vacuum Sublimation

Many organic compounds sublime readily under reduced pressure. When the vapor pressure of the solid equals the applied pressure, sublimation occurs, and the behavior is identical to that shown in Figure 17.1B. The solid phase passes directly into the vapor phase. From the data given in Table 17.1, you should expect camphor, naphthalene, and benzoic acid to sublime at or below the respective applied pressures of 370, 7, and 6 mm Hg. In principle, you can sublime *p*-nitrobenzaldehyde (last entry in the table), but it would not be practical because of the low applied pressure required.

17.4 Sublimation Methods

Sublimation can be used to purify solids. The solid is warmed until its vapor pressure becomes high enough for it to vaporize and condense as a solid on a cooled surface placed closely above. Several types of apparatus are illustrated in Figure 17.2. In each case, the cooled condensing surface is a tube filled with ice-cold water. The tube is filled from a beaker containing ice and water by using a Pasteur pipette. If the cooling water becomes warm before the sublimation is completed, the tube is emptied and refilled, once again by using a Pasteur pipette for these operations.

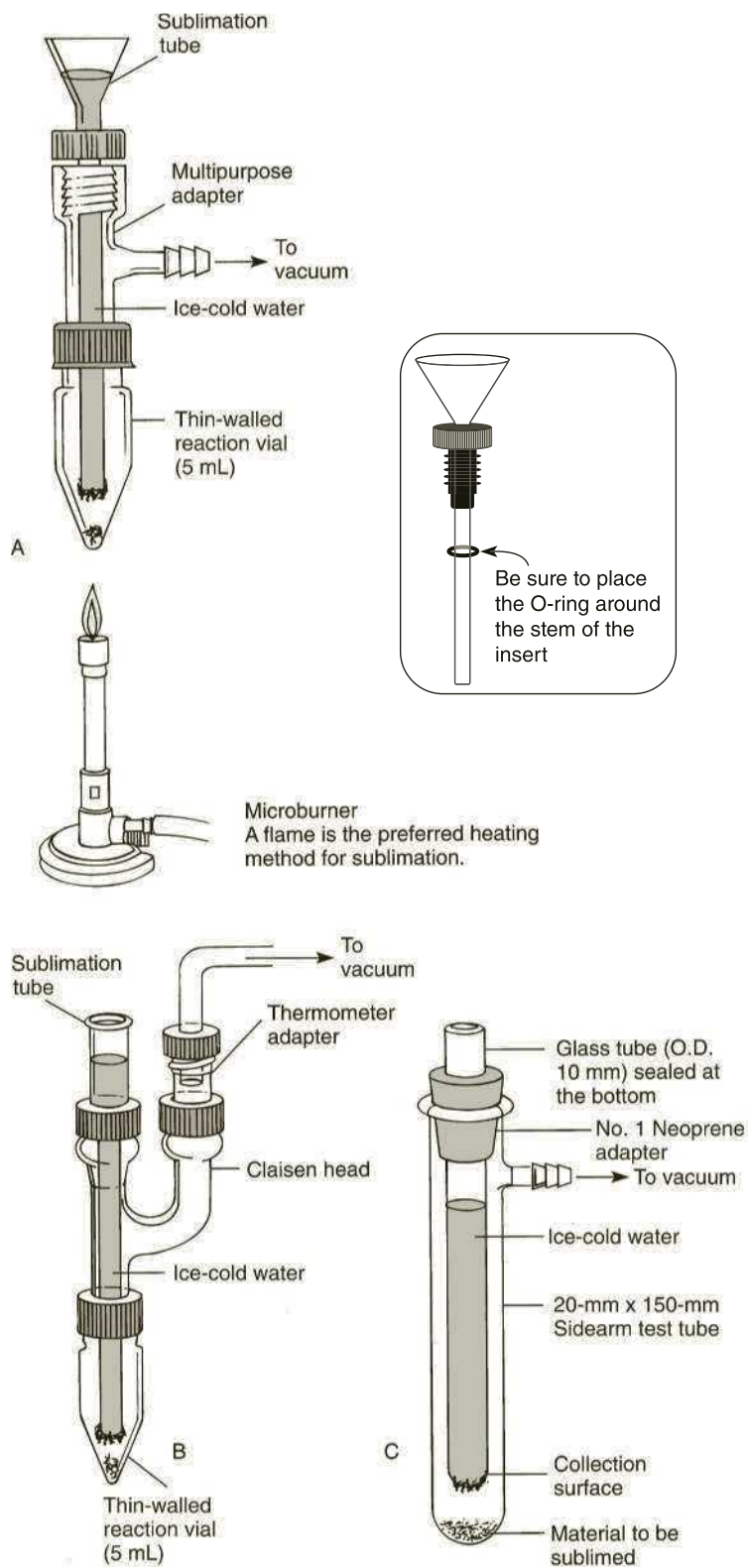


Figure 172
Sublimation apparatus.

A flame is the preferred heating device for a sublimation. The burner can be held by its cool base (not the hot barrel!) and moved up and down the sides of the thin-walled outer vial or tube to “chase” any solid that has formed on the sides toward the cold tube in the center. With an aluminum block, a ring of solid often forms on the inside walls of the apparatus just where it leaves the heating block. If this happens, using the aluminum collars will improve the situation considerably. When using a conical vial, use a thin-walled conical vial instead of a regular conical vial, because the thicker glass can shatter when heated by a flame.

Many solids do not develop enough vapor pressure at atmospheric pressure (760 mm Hg) to be purified by sublimation, but they frequently can be sublimed at reduced pressure. Thus, most sublimation equipment has provision for connection to an aspirator or other vacuum source. Reduction of pressure also helps to prevent thermal decomposition of substances that would require high temperatures to sublime at ordinary pressures.

Remember that while performing a sublimation, it is important to keep the temperature below the melting point of the solid. After sublimation, the material that has collected on the cooled surface is recovered by removing the central tube (cold-finger) from the apparatus. Take care in removing this tube to avoid dislodging the crystals that have collected. The deposit of crystals is scraped from the inner tube with a spatula. If reduced pressure has been used, the pressure must be released carefully to keep a blast of air from dislodging the crystals.

17.5 Advantages of Sublimation

One advantage of sublimation is that no solvent is used, and therefore none needs to be removed later. Sublimation also removes occluded material, such as molecules of solvent, from the sublimed substance. For instance, caffeine (sublimes at 178°C, melts at 236°C) absorbs water gradually from the atmosphere to form a hydrate. During sublimation, this water is lost, and anhydrous caffeine is obtained. If too much solvent is present in a sample to be sublimed, however, instead of becoming lost, it condenses on the cooled surface and thus interferes with the sublimation.

Sublimation is a faster method of purification than crystallization but not as selective. Similar vapor pressures are often a factor in dealing with solids that sublime; consequently, little separation can be achieved. For this reason, solids are far more often purified by crystallization. Sublimation is most effective in removing a volatile substance from a nonvolatile compound, particularly a salt or other inorganic material. Sublimation is also effective in removing highly volatile bicyclic or other symmetrical molecules from less volatile reaction products. Examples of volatile bicyclic compounds are borneol, isoborneol, and camphor.

PROBLEMS

1. Why is solid carbon dioxide called dry ice? How does it differ from solid water in behavior?
2. Under what conditions can you have *liquid* carbon dioxide?
3. A solid substance has a vapor pressure of 800 mm Hg at its melting point (80°C). Describe how the solid behaves as the temperature is raised from room temperature to 80°C while the atmospheric pressure is held constant at 760 mm Hg.
4. A solid substance has a vapor pressure of 100 mm Hg at the melting point (100°C). Assuming an atmospheric pressure of 760 mm Hg, describe the behavior of this solid as the temperature is raised from room temperature to its melting point.
5. A substance has a vapor pressure of 50 mm Hg at the melting point (100°C). Describe how you would experimentally sublime this substance.

Steam Distillation

The simple, vacuum, and fractional distillations described in Techniques 14, 15, and 16 are applicable to completely soluble (miscible) mixtures only. When liquids are *not* mutually soluble (immiscible), they can also be distilled, but with a somewhat different result. A mixture of immiscible liquids will boil at a lower temperature than the boiling points of any of the separate components as pure compounds. When steam is used to provide one of the immiscible phases, the process is called **steam distillation**. The advantage of this technique is that the desired material distills at a temperature below 100°C. Thus, if unstable or very high-boiling substances are to be removed from a mixture, decomposition is avoided. Because all gases mix, the two substances can mix in the vapor and codistill. Once the distillate is cooled, the desired component, which is not miscible, separates from the water. Steam distillation is used widely in isolating liquids from natural sources. It is also used in removing a reaction product from a tarry reaction mixture.

18.1 Differences between Distillation of Miscible and Immiscible Mixtures

Two liquids, A and B, that are mutually soluble (miscible), and that do not interact, form an ideal solution and follow Raoult's Law, as shown in Equation 1.

$$\text{Miscible Liquids} \quad P_{\text{total}} = P_A^\circ N_A + P_B^\circ N_B \quad [1]$$


Note that the vapor pressures of pure liquids P_A° and P_B° are not added directly to give the total pressure P_{total} but are reduced by the respective mole fractions N_A and N_B . The total pressure above a miscible or homogeneous solution will depend on P_A° and P_B° and also N_A and N_B . Thus, the composition of the vapor will also depend on *both* the vapor pressures and the mole fractions of each component.

$$\text{Immiscible Liquids} \quad P_{\text{total}} = P_A^\circ + P_B^\circ \quad [2]$$

In contrast, when two mutually insoluble (immiscible) liquids are "mixed" to give a heterogeneous mixture, each exerts its own vapor pressure, independently of the other, as shown in Equation 2. The mole fraction term does not appear in this equation, because the compounds are not miscible. You simply add the vapor pressures of the pure liquids P_A° and P_B° at a given temperature to obtain the total pressure above the mixture. When the total pressure equals 760 mm Hg, the mixture boils. The composition of the vapor from an immiscible mixture, in contrast to the miscible mixture, is determined only by the vapor pressures of the two substances codistilling. Equation 3 defines the composition of the vapor from an immiscible mixture. Calculations involving this equation are given in Section 18.2.

$$\frac{\text{Moles A}}{\text{Moles B}} = \frac{P_A^\circ}{P_B^\circ} \quad [3]$$

A mixture of two immiscible liquids boils at a lower temperature than the boiling points of either component. The explanation for this behavior is like that given

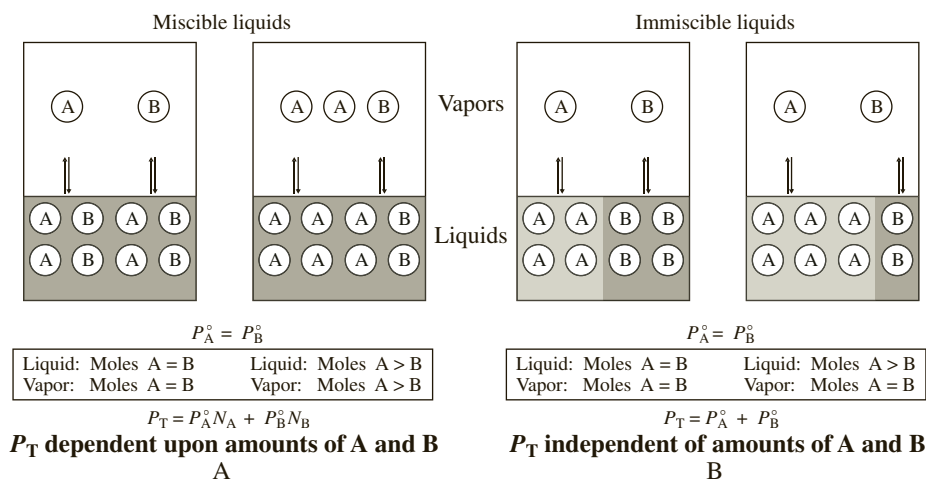**Figure 18.1**

Total pressure behavior for miscible and immiscible liquids. (A) Ideal miscible liquids follow Raoult's Law: P_T depends on the mole fractions and vapor pressures of A and of B. (B) Immiscible liquids do not follow Raoult's Law: P_T depends only on the vapor pressures of A and B.

for minimum-boiling-point azeotropes (Technique 15, Section 15.7). Immiscible liquids behave as they do because an extreme incompatibility between the two liquids leads to higher combined vapor pressures than Raoult's Law would predict. The higher combined vapor pressures cause a lower boiling point for the mixture than for either single component. Thus, you may think of steam distillation as a special type of azeotropic distillation in which the substance is completely insoluble in water.

The differences in behavior of miscible and immiscible liquids, where it is assumed that P_A° equals P_B° , are shown in Figure 18.1. Note that with miscible liquids, the composition of the vapor depends on the relative amounts of A and B present (Figure 18.1A). Thus, the composition of the vapor must change during a distillation. In contrast, the composition of the vapor with immiscible liquids is independent of the amounts of A and B present (Figure 18.1B). Hence, the vapor composition must remain *constant* during the distillation of such liquids, as predicted by Equation 3. Immiscible liquids act as if they were being distilled simultaneously from separate compartments, as shown in Figure 18.1B, even though in practice they are "mixed" during a steam distillation. Because all gases mix, they do give rise to a homogeneous vapor and codistill.

18.2 Immiscible Mixtures: Calculations

The composition of the distillate is constant during a steam distillation, as is the boiling point of the mixture. The boiling points of steam-distilled mixtures will always be below the boiling point of water (bp 100°C), as well as the boiling point of any of the other substances distilled. Some representative boiling points and compositions of steam distillates are given in Table 18.1. Note that the higher the boiling point of a pure substance, the more closely the temperature of the steam distillate approaches, but does not exceed, 100°C. This is a reasonably low temperature, and it avoids the decomposition that might result at high temperatures with a simple distillation.

For immiscible liquids, the molar proportions of two components in a distillate equal the ratio of their vapor pressures in the boiling mixture, as given in Equation 3. When Equation 3 is rewritten for an immiscible mixture involving water,

Equation 4 results. Equation 4 can be modified by substituting the relation moles = (weight/molecular weight) to give Equation 5.

$$\frac{\text{Moles substance}}{\text{Moles water}} = \frac{P_{\text{substance}}^{\circ}}{P_{\text{water}}^{\circ}} \quad [4]$$

$$\frac{\text{Wt substance}}{\text{Wt water}} = \frac{(P_{\text{substance}}^{\circ})(\text{molecular weight}_{\text{substance}})}{(P_{\text{water}}^{\circ})(\text{molecular weight}_{\text{water}})} \quad [5]$$

A sample calculation using this equation is given in Figure 18.2. Notice that the result of this calculation is very close to the experimental value given in Table 18.1.

18.3 Steam Distillation: Methods

Two methods for steam distillation are in general use in the laboratory: the **direct method** and the **live steam method**. In the first method, steam is generated *in situ* (in place) by heating a distillation flask containing the compound and water. In the second method, steam is generated outside and is passed into the distillation flask using an inlet tube.

Problem How many grams of water must be distilled to steam-distill 1.55 g of 1-octanol from an aqueous solution? What will be the composition (wt %) of the distillate? The mixture distills at 99.4°C.

Answer The vapor pressure of water at 99.4°C must be obtained from the *CRC Handbook* (= 744 mm Hg).

(a) Obtain the partial pressure of 1-octanol.

$$P_{\text{1-octanol}}^{\circ} = P_{\text{total}} - P_{\text{water}}^{\circ}$$

$$P_{\text{1-octanol}}^{\circ} = (760 - 744) = 16 \text{ mm Hg}$$

(b) Obtain the composition of the distillate.

$$\frac{\text{wt 1-octanol}}{\text{wt water}} = \frac{(16)(130)}{(744)(18)} = 0.155 \text{ g/g-water}$$

(c) Clearly, 10 g of water must be distilled.

$$(0.155 \text{ g/g-water})(10 \text{ g-water}) = 1.55 \text{ g 1-octanol}$$

(d) Calculate the **weight** percentages.

$$\text{1-octanol} = 1.55 \text{ g}/(10 \text{ g} + 1.55 \text{ g}) = 13.4\%$$

$$\text{water} = 10 \text{ g}/(10 \text{ g} + 1.55 \text{ g}) = 86.6\%$$

Figure 18.2

Sample calculations for a steam distillation.

Table 18.1 Boiling points and compositions of steam distillates

Mixture	Boiling Point of Pure Substance (°C)	Boiling Point of Mixture (°C)	Composition (% Water)
Benzene–water	80.1	69.4	8.9%
Toluene–water	110.6	85.0	20.2%
Hexane–water	69.0	61.6	5.6%
Heptane–water	98.4	79.2	12.9%
Octane–water	125.7	89.6	25.5%
Nonane–water	150.8	95.0	39.8%
1-Octanol–water	195.0	99.4	90.0%

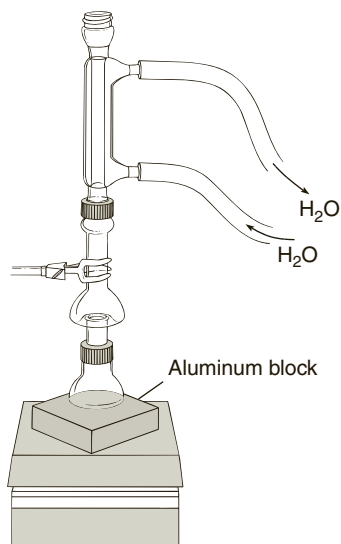


Figure 18.3
Microscale steam distillation.

A. Direct Method

Microscale. The direct method of steam distillation is the only one suitable for microscale reactions. Steam is produced in the conical vial or distillation flask (*in situ*) by heating water to its boiling point in the presence of the compound to be distilled. This method works well for small amounts of materials. A microscale steam distillation apparatus is shown in Figure 18.3. Water and the compound to be distilled are placed in the flask and heated. A stirring bar or a boiling stone should be used to prevent bumping. The vapors of the water and the desired compound codistill when they are heated. They are condensed and collect in the Hickman head. When the Hickman head fills, the distillate is removed with a Pasteur pipette and placed in another vial for storage. For the typical microscale experiment, it will be necessary to fill the well and remove the distillate three or four times. All these distillate fractions are placed in the same storage container. The efficiency in collecting the distillate can sometimes be improved if the inside walls of the Hickman head are rinsed several times into the well. A Pasteur pipette is used to perform the rinsing. Distillate is withdrawn from the well, and then it is used to wash the walls of the Hickman head all the way around the head. After the walls have been washed and when the well is full, the distillate can be withdrawn and transferred to the storage container. It may be necessary to add more water during the course of the distillation. More water is added (remove the condenser if used) through the center of the Hickman head by using a Pasteur pipette.

Semimicroscale. The apparatus shown in Technique 14, Figure 14.10, may also be used to perform a steam distillation at the microscale level or slightly above. This apparatus avoids the need to empty the collected distillate during the course of the distillation as is required when a Hickman head is used.

Macroscale. A larger-scale direct method steam distillation is illustrated in Figure 18.4. Although a heating mantle may be used, it is probably best to use a flame with this

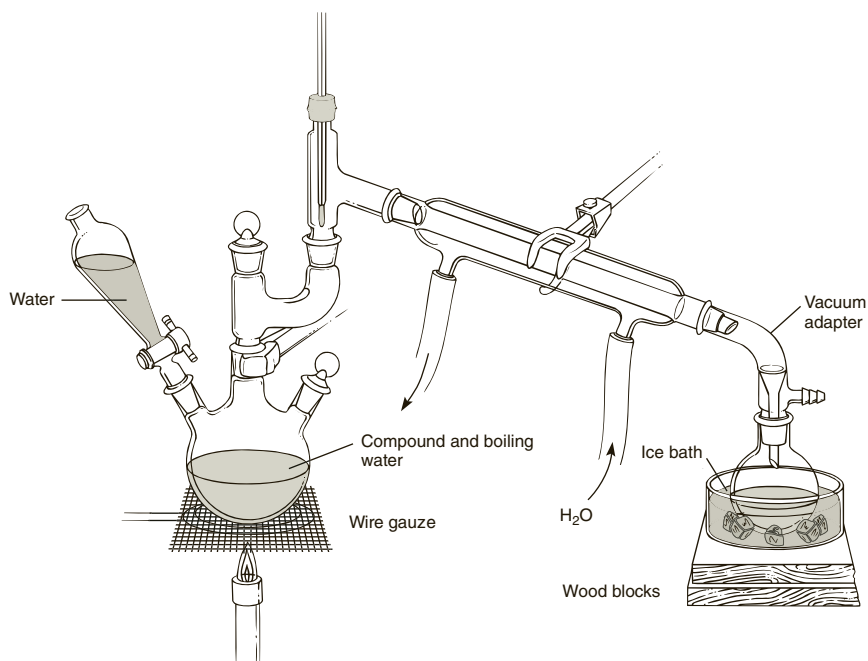


Figure 18.4
Macroscale direct steam distillation.

method because a large volume of water must be heated rapidly. A boiling stone must be used to prevent bumping. The separatory funnel allows more water to be added during the course of the distillation.

B. Live Steam Method

Macroscale. A large-scale steam distillation using the live steam method is shown in Figure 18.5. If steam lines are available in the laboratory, they may be attached directly to the steam trap (purge them first to drain water). If steam lines are not available, an external steam generator (see inset) must be prepared. The external generator usually will require a flame to produce steam at a rate fast enough for the distillation. When the distillation is first started, the clamp at the bottom of the steam trap is left open. The steam lines will have a large quantity of condensed water in them until they are well heated. When the lines become hot and condensation of steam ceases, the clamp may be closed. Occasionally, the clamp will have to be reopened to remove condensate. In this method, the steam agitates the mixture as it enters the bottom of the flask, and a stirrer or boiling stone is not required.

CAUTION



Hot steam can produce severe burns.

Sometimes it is helpful to heat the three-necked distilling flask with a heating mantle (or flame) to prevent excessive condensation at that point. Steam must be

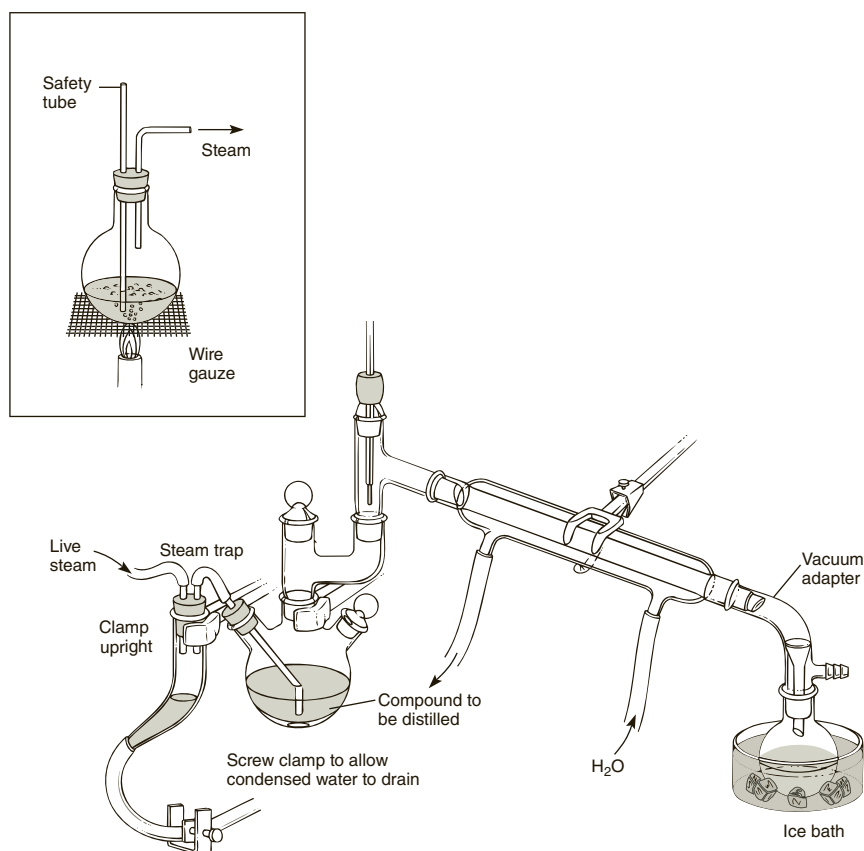


Figure 18.5
Macroscale steam distillation using live steam.

admitted fast enough for you to see the distillate condensing as a milky white fluid in the condenser. The vapors that codistill will separate on cooling to give this cloudiness. When the condensate becomes clear, the distillation is near the end. The flow of water through the condenser should be faster than in other types of distillation to help cool the vapors. Make sure the vacuum adapter remains cool to the touch. An ice bath may be used to cool the receiving flask if desired. When the distillation is to be stopped, the screw clamp on the steam trap should be opened, and the steam inlet tube must be removed from the three-necked flask. If this is not done, liquid will back up into the tube and steam trap.

PROBLEMS

1. Calculate the weight of benzene codistilled with each gram of water and the percentage composition of the vapor produced during a steam distillation. The boiling point of the mixture is 69.4°C. The vapor pressure of water at 69.4°C is 227.7 mm Hg. Compare the result with the data in Table 18.1.
2. Calculate the approximate boiling point of a mixture of bromobenzene and water at atmospheric pressure. A table of vapor pressures of water and bromobenzene at various temperatures is given.

Temperature (°C)	Vapor Pressures (mm Hg)	
	Water	Bromobenzene
93	588	110
94	611	114
95	634	118
96	657	122
97	682	127
98	707	131
99	733	136

3. Calculate the weight of nitrobenzene that codistills (bp 99°C) with each gram of water during a steam distillation. You may need the data given in Problem 2.
4. A mixture of *p*-nitrophenol and *o*-nitrophenol can be separated by steam distillation. The *o*-nitrophenol is steam volatile, and the *para* isomer is not volatile. Explain. Base your answer on the ability of the isomers to form hydrogen bonds internally.
5. When another compound is mixed with water and distilled, how can you determine if it is a steam distillation or an azeotropic distillation?

Column Chromatography


The most modern and sophisticated methods of separating mixtures that the organic chemist has available all involve **chromatography**. Chromatography is defined as the separation of a mixture of two or more compounds or ions by distribution between two phases, one of which is stationary and the other moving. Various types of chromatography are possible, depending on the nature of the two phases involved: **solid–liquid** (column, thin-layer, and paper), **liquid–liquid** (high-performance liquid), and **gas–liquid** (vapor-phase) chromatographic methods are common.

All chromatography works on much the same principle as solvent extraction (Technique 12). Basically, the methods depend on the differential solubilities or adsorptivities of the substances to be separated relative to the two phases between which they are to be partitioned. In this chapter, column chromatography, a solid–liquid method, is considered. Thin-layer chromatography is examined in Technique 20; high-performance liquid chromatography is discussed in Technique 21; and gas chromatography, a gas–liquid method, is discussed in Technique 22.

19.1 Adsorbents

Column chromatography is a technique based on both adsorptivity and solubility. It is a solid–liquid phase-partitioning technique. The solid may be almost any material that does not dissolve in the associated liquid phase; the solids used most commonly are silica gel, $\text{SiO}_2 \cdot x\text{H}_2\text{O}$, also called silicic acid, and alumina, $\text{Al}_2\text{O}_3 \cdot x\text{H}_2\text{O}$. These compounds are used in their powdered or finely ground forms (usually 200 to 400 mesh).¹

Most alumina used for chromatography is prepared from the impure ore bauxite $\text{Al}_2\text{O}_3 \cdot x\text{H}_2\text{O} + \text{Fe}_2\text{O}_3$. The bauxite is dissolved in hot sodium hydroxide and filtered to remove the insoluble iron oxides; the alumina in the ore forms the soluble amphoteric hydroxide $\text{Al}(\text{OH})_4^-$. The hydroxide is precipitated by CO_2 , which reduces the pH, as $\text{Al}(\text{OH})_3$. When heated, the $\text{Al}(\text{OH})_3$ loses water to form pure alumina Al_2O_3 .

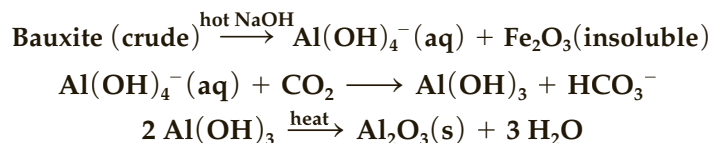

Alumina prepared in this way is called **basic alumina** because it still contains some hydroxides. Basic alumina cannot be used for chromatography of compounds that are base sensitive. Therefore, it is washed with acid to neutralize the base,

¹ The term *mesh* refers to the number of openings per linear inch found in a screen. A large number refers to a fine screen (finer wires more closely spaced). When particles are sieved through a series of these screens, they are classified by the smallest mesh screen that they will pass through. Mesh 5 would represent a coarse gravel, and mesh 800 would be a fine powder.

giving **acid-washed alumina**. This material is unsatisfactory unless it has been washed with enough water to remove *all* the acid; on being so washed, it becomes the best chromatographic material, called **neutral alumina**. If a compound is acid-sensitive, either basic or neutral alumina must be used. You should be careful to ascertain what type of alumina is being used for chromatography. Silica gel is not available in any form other than that suitable for chromatography.

19.2 Interactions

If powdered or finely ground alumina (or silica gel) is added to a solution containing an organic compound, some of the organic compound will **adsorb** onto or adhere to the fine particles of alumina. Many kinds of intermolecular forces cause organic molecules to bind to alumina. These forces vary in strength according to their type. Nonpolar compounds bind to the alumina using only van der Waals forces. These are weak forces, and nonpolar molecules do not bind strongly unless they have extremely high molecular weights. The most important interactions are those typical of polar organic compounds. Either these forces are of the dipole–dipole type or they involve some direct interaction (coordination, hydrogen bonding, or salt formation). These types of interactions are illustrated in Figure 19.1, which for convenience shows only a portion of the alumina structure. Similar interactions occur with silica gel. The strengths of such interactions vary in the approximate order:

Salt formation > coordination > hydrogen-bonding > dipole–dipole > van der Waals

Strength of interaction varies among compounds. For instance, a strongly basic amine would bind more strongly than a weakly basic one (by coordination). In fact, strong bases and strong acids often interact so strongly that they *dissolve* alumina to some extent. You can use the following rule of thumb:

NOTE: The more polar the functional group, the stronger the bond to alumina (or silica gel).

A similar rule holds for solubility. Polar solvents dissolve polar compounds more effectively than nonpolar solvents; nonpolar compounds are dissolved best by nonpolar solvents. Thus, the extent to which any given solvent can wash an adsorbed compound from alumina depends almost directly on the relative polarity of the solvent. For example, although a ketone adsorbed on alumina might not be

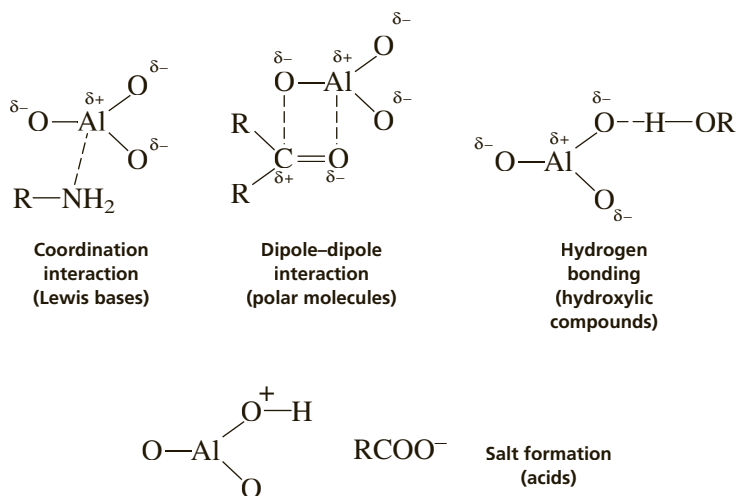


Figure 19.1
Possible interactions of organic compounds with alumina.

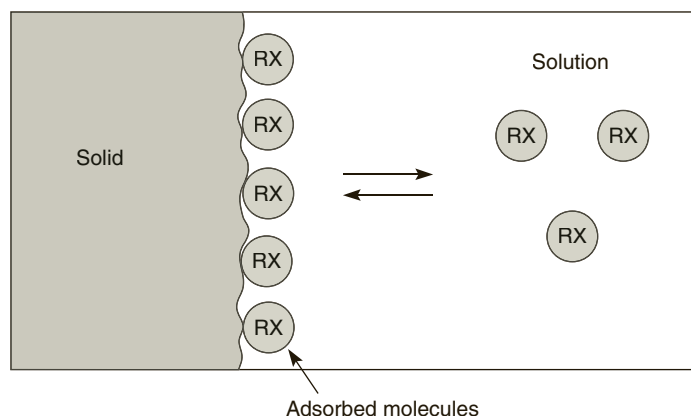


Figure 19.2
Dynamic adsorption equilibrium.

removed by hexane, it might be removed completely by chloroform. For any adsorbed material, a kind of **distribution** equilibrium can be envisioned between the adsorbent material and the solvent. This is illustrated in Figure 19.2.

The distribution equilibrium is *dynamic*, with molecules constantly *adsorbing* from the solution and *desorbing* into it. The average number of molecules remaining adsorbed on the solid particles at equilibrium depends both on the particular molecule (RX) involved and the dissolving power of the solvent with which the adsorbent must compete.

19.3 Principle of Column Chromatographic Separation

The dynamic equilibrium mentioned previously, and the variations in the extent to which different compounds adsorb on alumina or silica gel, underlie a versatile and ingenious method for *separating* mixtures of organic compounds. In this method, the mixture of compounds to be separated is introduced onto the top of a cylindrical glass column (Figure 19.3) *packed*, or filled, with fine alumina particles (stationary solid phase). The adsorbent is continuously washed by a flow of solvent (moving phase) passing through the column.

Initially, the components of the mixture adsorb onto the alumina particles at the top of the column. The continuous flow of solvent through the column **elutes**, or washes, the solutes off the alumina and sweeps them down the column. The solutes (or materials to be separated) are called **eluates** or **elutants**, and the solvents are called **eluent**s. As the solutes pass down the column to fresh alumina, new equilibria are established among the adsorbent, the solutes, and the solvent. The constant equilibration means that different compounds will move down the column at differing rates, depending on their relative affinity for the adsorbent on one hand and for the solvent on the other. Because the number of alumina particles is large, because they are closely packed, and because fresh solvent is being added continuously, the number of equilibrations between adsorbent and solvent that the solutes experience is enormous.

As the components of the mixture are separated, they begin to form moving bands (or zones), each band containing a single component. If the column is long enough and the other parameters (column diameter, adsorbent, solvent, and flow rate) are correctly chosen, the bands separate from one another, leaving gaps of pure solvent in between. As each band (solvent and solute) passes out the bottom of the column, it can be collected before the next band arrives. If the parameters mentioned are poorly chosen, the

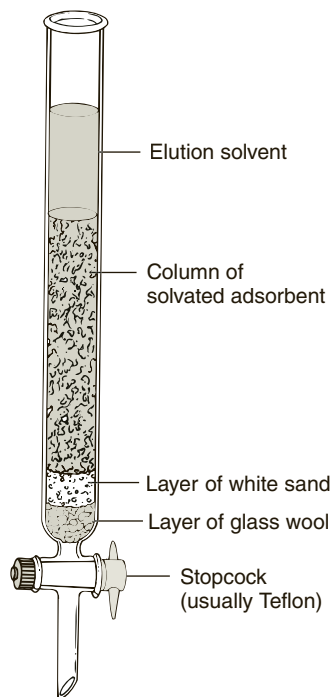


Figure 19.3
Chromatographic column.

various bands either overlap or coincide, in which case either a poor separation or no separation at all is the result. A successful chromatographic separation is illustrated in Figure 19.4.

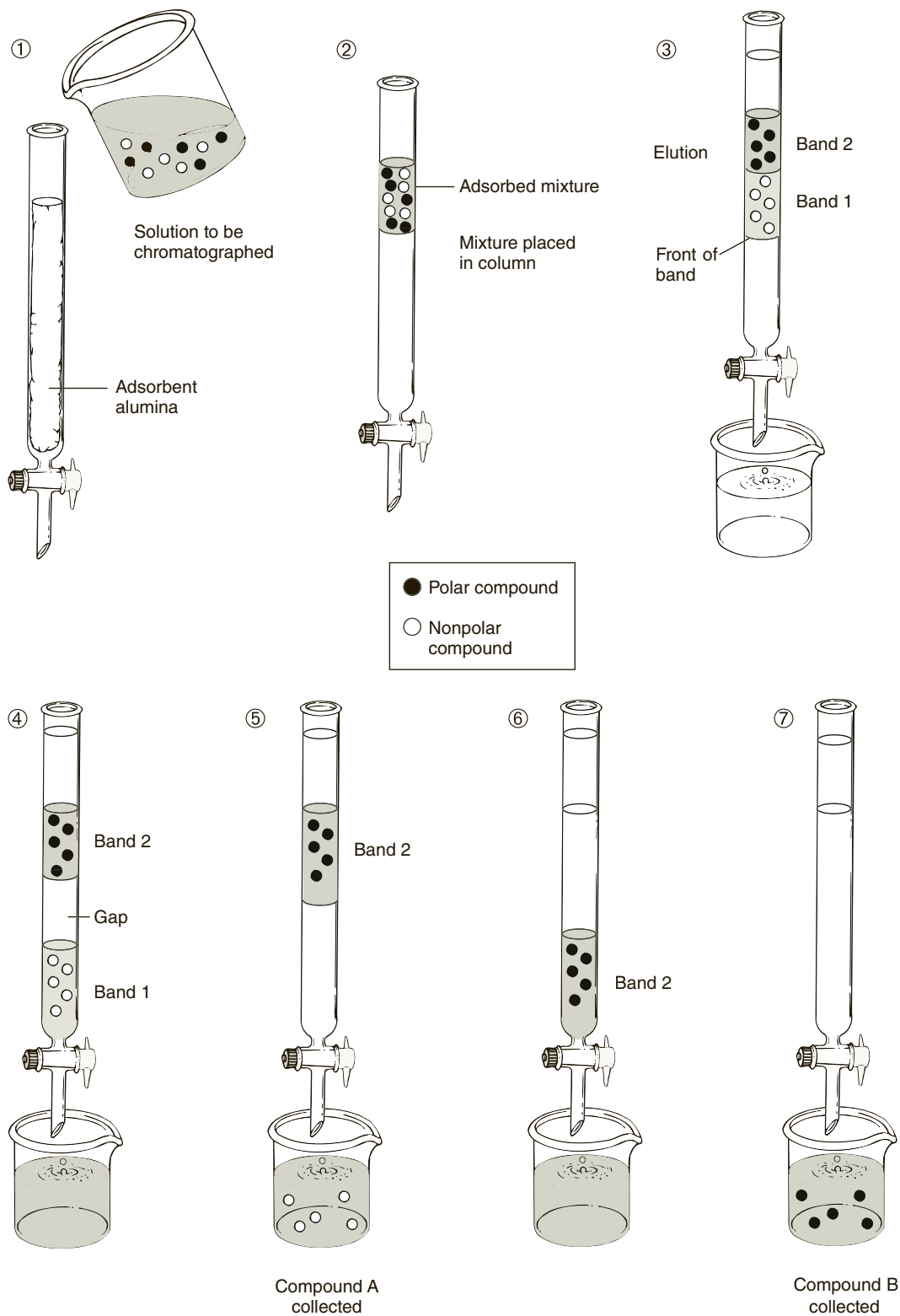


Figure 19.4
Sequence of steps in a chromatographic separation.

19.4 Parameters Affecting Separation

The versatility of column chromatography results from the many factors that can be adjusted. These include

1. Adsorbent chosen
2. Polarity of the solvents chosen
3. Size of the column (both length and diameter) relative to the amount of material to be chromatographed
4. Rate of elution (or flow)

By careful choosing of the conditions, almost any mixture can be separated. This technique has even been used to separate optical isomers. An optically active solid-phase adsorbent was used to separate the enantiomers.

Two fundamental choices for anyone attempting a chromatographic separation are the kind of adsorbent and the solvent system. In general, nonpolar compounds pass through the column faster than polar compounds because they have a smaller affinity for the adsorbent. If the adsorbent chosen binds all the solute molecules (both polar and nonpolar) strongly, they will not move down the column. On the other hand, if too polar a solvent is chosen, all the solutes (polar and nonpolar) may simply be washed through the column, with no separation taking place. The adsorbent and the solvent should be chosen so that neither is favored excessively in the equilibrium competition for solute molecules.²

A. Adsorbents

In Table 19.1, various kinds of adsorbents (solid phases) used in column chromatography are listed. The choice of adsorbent often depends on the types of compounds to be separated. Cellulose, starch, and sugars are used for polyfunctional plant and animal materials (natural products) that are very sensitive to acid–base interactions. Magnesium silicate is often used for separating acetylated sugars, steroids, and essential oils. Silica gel and Florisil are relatively mild toward most compounds and

Table 19.1 Solid adsorbents for column chromatography

Paper	Increasing strength of binding interactions toward polar compounds ↓
Cellulose	
Starch	
Sugars	
Magnesium silicate	
Calcium sulfate	
Silicic acid	
Silica gel	
Florisil	
Magnesium oxide	
Aluminum oxide (alumina)*	
Activated charcoal (Norit)	

*Basic, acid-washed, and neutral.

² Often the chemist uses thin-layer chromatography (TLC), which is described in Technique 20, to arrive at the best choices of solvents and adsorbents for the best separation. The TLC experimentation can be performed quickly and with extremely small amounts (microgram quantities) of the mixture to be separated. This saves significant time and materials. Technique 20 describes this use of TLC.

are widely used for a variety of functional groups—hydrocarbons, alcohols, ketones, esters, acids, azo compounds, and amines. Alumina is the most widely used adsorbent and is obtained in the three forms mentioned in Section 19.1: acidic, basic, and neutral. The pH of acidic or acid-washed alumina is approximately 4. This adsorbent is particularly useful for separating acidic materials such as carboxylic acids and amino acids. Basic alumina has a pH of 10 and is useful in separating amines. Neutral alumina can be used to separate a variety of nonacidic and nonbasic materials.


The approximate strength of the various adsorbents listed in Table 19.1 is also given. The order is only approximate, and therefore it may vary. For instance, the strength, or separating abilities, of alumina and silica gel largely depend on the amount of water present. Water binds tightly to either adsorbent, taking up sites on the particles that could otherwise be used for equilibration with solute molecules. If one adds water to the adsorbent, it is said to have been **deactivated**. Anhydrous alumina or silica gel is said to be highly **activated**. High activity is usually avoided with these adsorbents. Use of the highly active forms of either alumina or silica gel, or of the acidic or basic forms of alumina, can often lead to molecular rearrangement or decomposition in certain types of solute compounds.

The chemist can select the degree of activity that is appropriate to carry out a particular separation. To accomplish this, highly activated alumina is mixed thoroughly with a precisely measured quantity of water. The water partially hydrates the alumina and thus reduces its activity. By carefully determining the amount of water required, the chemist can have available an entire spectrum of possible activities.

B. Solvents

In Table 19.2, some common chromatographic solvents are listed, along with their relative ability to dissolve polar compounds. Sometimes a single solvent can be found that will separate all the components of a mixture. Sometimes a mixture of solvents can be found that will achieve separation. More often, you must start elution with a nonpolar solvent to remove relatively nonpolar compounds from the column and then gradually increase the solvent polarity to force compounds of greater polarity to come down the column, or elute. The approximate order in

Table 19.2 Solvents (eluent) for chromatography

Petroleum ether	Increasing polarity and “solvent power” toward polar functional groups 
Cyclohexane	
Carbon tetrachloride*	
Toluene	
Chloroform*	
Methylene chloride	
Diethyl ether	
Ethyl acetate	
Acetone	
Pyridine	
Ethanol	
Methanol	
Water	
Acetic acid	

*Suspected carcinogens.


Table 19.3 Elution sequence for compounds

Hydrocarbons	Fastest (will elute with nonpolar solvent) Order of elution Slowest (need a polar solvent)
Olefins	
Ethers	
Halocarbons	
Aromatics	
Ketones	
Aldehydes	
Esters	
Alcohols	
Amines	
Acids, strong bases	▼

which various classes of compounds elute by this procedure is given in Table 19.3. In general, nonpolar compounds travel through the column faster (elute first), and polar compounds travel more slowly (elute last). However, molecular weight is also a factor in determining the order of elution. A nonpolar compound of high molecular weight travels more slowly than a nonpolar compound of low molecular weight, and it may even be passed by some polar compounds.

Solvent polarity functions in two ways in column chromatography. First, a polar solvent will better dissolve a polar compound and move it down the column faster. Therefore, as already mentioned, one usually increases the polarity of the solvent during column chromatography to wash down compounds of increasing polarity. Second, as the polarity of the solvent increases, the solvent itself will displace adsorbed molecules from the alumina or silica and take their place on the column. Because of this second effect, a polar solvent will move *all* types of compounds, both polar and nonpolar, down the column at a faster rate than a nonpolar solvent.

When the polarity of the solvent has to be changed during a chromatographic separation, some precautions must be taken. Rapid changes from one solvent to another are to be avoided (especially when silica gel or alumina is involved). Usually, small percentages of a new solvent are mixed slowly into the one in use until the percentage reaches the desired level. If this is not done, the column packing often “cracks” as a result of the heat liberated when alumina or silica gel is mixed with a solvent. The solvent solvates the adsorbent, and the formation of a weak bond generates heat.

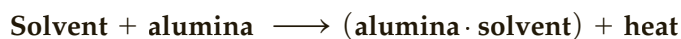

Often, enough heat is generated locally to evaporate the solvent. The formation of vapor creates bubbles, which forces a separation of the column packing; this is called **cracking**. A cracked column does not produce a good separation, because it has discontinuities in the packing. The way in which a column is packed or filled is also very important in preventing cracking.

Certain solvents should be avoided with alumina or silica gel, especially with the acidic, basic, and highly active forms. For instance, with any of these adsorbents, acetone dimerizes via an aldol condensation to give diacetone alcohol. Mixtures of esters *transesterify* (exchange their alcoholic portions) when ethyl acetate or an alcohol is the eluent. Finally, the most active solvents (pyridine, methanol, water, and acetic acid) dissolve and elute some of the adsorbent itself. Generally, try to avoid going to solvents more polar than diethyl ether or methylene chloride in the eluent series (Table 19.2).

Table 19.4 Size of column and amount of adsorbent for typical sample sizes

Amount of Sample (g)	Amount of Adsorbent (g)	Column Diameter (mm)	Column Height (mm)
0.01	0.3	3.5	30
0.10	3.0	7.5	60
1.00	30.0	16.0	130
10.00	300.0	35.0	280

C. Column Size and Adsorbent Quantity

The column size and the amount of adsorbent must also be selected correctly to separate a given amount of sample well. As a rule of thumb, the amount of adsorbent should be 25 to 30 times, by weight, the amount of material to be separated by chromatography. Furthermore, the column should have a height-to-diameter ratio of about 8:1. Some typical relations of this sort are given in Table 19.4.

Note, as a caution, that the difficulty of the separation is also a factor in determining the size and length of the column to be used and the amount of adsorbent needed. Compounds that do not separate easily may require longer columns and more adsorbent than specified in Table 19.4. For easily separated compounds, a shorter column and less adsorbent may suffice.

D. Flow Rate

The rate at which solvent flows through the column is also significant in the effectiveness of a separation. In general, the time the mixture to be separated remains on the column is directly proportional to the extent of equilibration between stationary and moving phases. Thus, similar compounds eventually separate if they remain on the column long enough. The time a material remains on the column depends on the flow rate of the solvent. If the flow is too slow, however, the dissolved substances in the mixture may diffuse faster than the rate at which they move down the column. Then the bands grow wider and more diffuse, and the separation becomes poor.

19.5 Packing the Column: Typical Problems

The most critical operation in column chromatography is packing (filling) the column with adsorbent. The **column packing** must be evenly packed and free of irregularities, air bubbles, and gaps. As a compound travels down the column, it moves in an advancing zone, or **band**. It is important that the leading edge, or **front**, of this band be horizontal, or perpendicular to the long axis of the column. If two bands are close together and do not have horizontal band fronts, it is impossible to collect one band while completely excluding the other. The leading edge of the second band begins to elute before the first band has finished eluting. This condition can be seen in Figure 19.5. There are two main reasons for this problem. First, if the top surface edge of the adsorbent packing is not level, nonhorizontal bands result. Second, bands may be nonhorizontal if the column is not held in an exactly vertical position in both planes (front to back and side to side). When you are preparing a column, you must watch both these factors carefully.

Another phenomenon, called **streaming** or **channeling**, occurs when part of the band front advances ahead of the major part of the band. Channeling occurs if there are any cracks or irregularities in the adsorbent surface or any irregularities caused by air bubbles in the packing. A part of the advancing front moves ahead of the rest of the band by flowing through the channel. Two examples of channeling are shown in Figure 19.6.

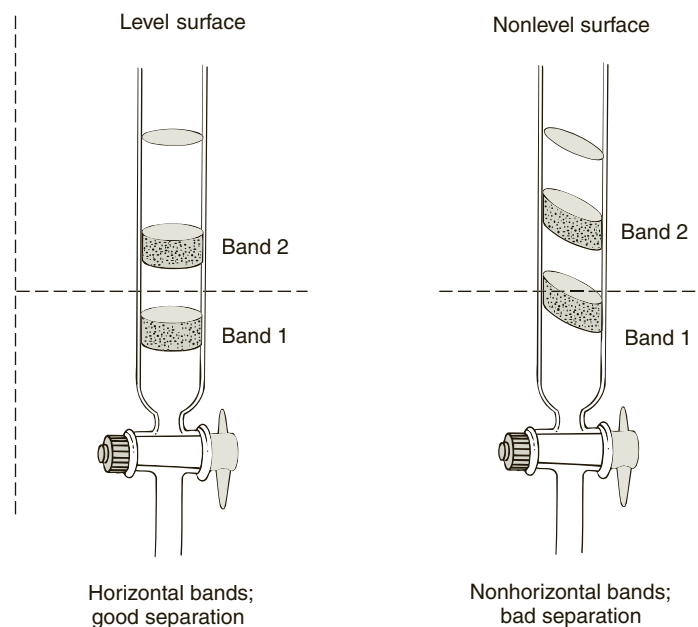


Figure 19.5
Comparison of horizontal and nonhorizontal band fronts.

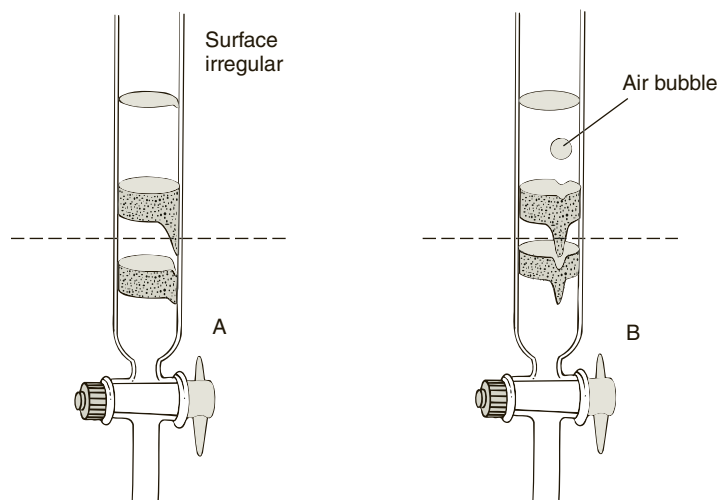


Figure 19.6
Channeling complications.

19.6 Packing the Column: Microscale Methods

The following methods are used to avoid problems resulting from uneven packing and column irregularities. These procedures should be followed carefully in preparing a chromatography column. Failure to pay close attention to the preparation of the column may well affect the quality of the separation.

Preparation of a column involves two distinct stages. In the first stage, a support base on which the packing will rest is prepared. This must be done so that the

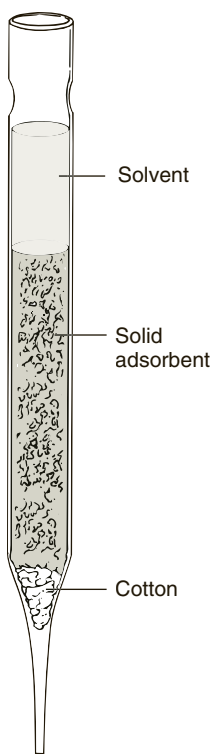


Figure 19.7
Microscale
chromatography
column.

packing, a finely divided material, does not wash out of the bottom of the column. In the second stage, the column of adsorbent is deposited on top of the supporting base.

A. Preparing the Support Base

For microscale applications, select a Pasteur pipette (5¾-inch) and clamp it upright (vertically). To reduce the amount of solvent needed to fill the column, break off most of the tip of the pipette. Place a small ball of cotton in the pipette and tamp it into position using a glass rod or a piece of wire. Take care not to plug the column totally by tamping the cotton too hard. The correct position of the cotton is shown in Figure 19.7. A microscale chromatography column is packed by one of the dry pack methods described in Part B of this section.

B. Depositing the Adsorbent

Dry Pack Method 1. To fill a microscale column, fill the Pasteur pipette (with the cotton plug, prepared as described in Section A) about half full with solvent. Using a microspatula, add the solid adsorbent slowly to the solvent in the column. As you add the solid, tap the column *gently* with a pencil, a finger, or a glass rod. The tapping promotes even settling and mixing and gives an evenly packed column free of air bubbles. As the adsorbent is added, solvent flows out of the Pasteur pipette. Because the adsorbent must not be allowed to dry during the packing process, you must use a means of controlling the solvent flow. If a piece of small-diameter plastic tubing is available, it can be fitted over the narrow tip of the Pasteur pipette. The flow rate can then be controlled using a screw clamp. A simple approach to controlling the flow rate is to use a finger over the top of the Pasteur pipette, much as you control the flow of liquid in a volumetric pipette. Continue adding the adsorbent slowly, with constant tapping, until the level of the adsorbent has reached the desired level. As you pack the column, be careful not to let the column run dry. The final column should appear as shown in Figure 19.7.

Dry Pack Method 2. An alternative dry pack method for microscale columns is to fill the Pasteur pipette with *dry* adsorbent, without any solvent. Position a plug of cotton in the bottom of the Pasteur pipette. The desired amount of adsorbent is added slowly, and the pipette tapped constantly, until the level of adsorbent has reached the desired height. Figure 19.7 can be used as a guide to judge the correct height of the column of adsorbent. When the column is packed, added solvent is allowed to percolate through the adsorbent until the entire column is moistened. The solvent is not added until just before the column is to be used.

NOTE: This method is not recommended for use with silica gel or for experiments where a careful separation is required.

This method is useful when the adsorbent is alumina, but it does not produce satisfactory results with silica gel. Even with alumina, poor separations can arise due to uneven packing, air bubbles, and cracking, especially if a solvent that has a highly exothermic heat of solvation is used.

19.7 Packing the Column: Semi-microscale and Macroscale Methods

As with microscale columns, the procedures described in this section should be followed carefully in preparing a semimicroscale or conventional-scale chromatography column. Failure to pay close attention to the details of these procedures may adversely affect the quality of the separation.

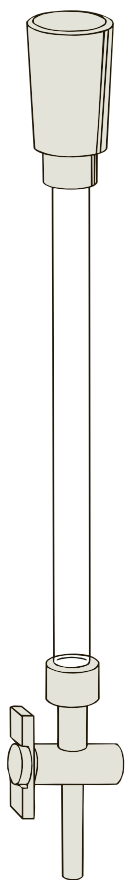


Figure 19.8
Commercial semimicroscale chromatography column. (The column is shown equipped with an optional solvent reservoir.)

Again, preparation of a column involves two distinct stages: preparing the support base and filling the column with adsorbent.

A. Preparing the Support Base

Semimicroscale Columns. An alternative apparatus for small-scale column chromatography is a commercial column, such as the one shown in Figure 19.8. This type of column is made of glass and has a solvent-resistant plastic stopcock at the bottom.³ The stopcock assembly contains a filter disc to support the adsorbent column. An optional upper fitting, also made of solvent-resistant plastic, serves as a solvent reservoir. The column shown in Figure 19.8 is equipped with the solvent reservoir. This type of column is available in a variety of lengths, ranging from 100 to 300 mm. Because the column has a built-in filter disc, it is not necessary to prepare a support base before the adsorbent is added.

Macroscopic Columns. For large-scale applications, clamp a chromatography column upright (vertically). The column (Figure 19.3) is a piece of cylindrical glass tubing with a stopcock attached at one end. The stopcock usually has a Teflon plug because stopcock grease (used on glass plugs) dissolves in many of the organic solvents used as eluents. Stopcock grease in the eluent will contaminate the eluates.

Instead of a stopcock, attach a piece of flexible tubing to the bottom of the column, with a screw clamp used to stop or regulate the flow (Figure 19.9). When a screw clamp is used, care must be taken that the tubing used is not dissolved by the solvents that will pass through the column during the experiment. Rubber, for instance, dissolves in chloroform, benzene, methylene chloride, toluene, or tetrahydrofuran (THF). Tygon tubing dissolves (actually, the plasticizer is removed) in many solvents, including benzene, methylene chloride, chloroform, ether, ethyl acetate, toluene, and THF. Polyethylene tubing is the best choice for use at the end of a column because it is inert with most solvents.

Next, the column is partially filled with a quantity of solvent, usually a nonpolar solvent like hexane, and a support for the finely divided adsorbent is prepared in the following way. A loose plug of glass wool is tamped down into the bottom of the column with a long glass rod until all entrapped air is forced out as bubbles. Take care not to plug the column totally by tamping the glass wool too hard. A small layer of clean white sand is formed on top of the glass wool by pouring sand into the column. The column is tapped to level the surface of the sand. Any sand adhering to the side of the column is washed down with a small quantity of solvent. The sand forms a base that supports the column of adsorbent and prevents it from washing through the stopcock. The column is packed in one of two ways: by the slurry method or by the dry pack method.

B. Depositing the Adsorbent

Slurry Method. The slurry method is not recommended as a microscale method for use with Pasteur pipettes. On a very small scale, it is too difficult to pack the column with the slurry without losing the solvent before the packing has been completed. Microscale columns should be packed by the dry pack method, as described in Section 19.6.

In the slurry method, the adsorbent is packed into the column as a mixture of a solvent and an undissolved solid. The slurry is prepared in a separate container

³ *Note to the Instructor:* With certain organic solvents, we have found that the “solvent-resistant” plastic stopcock may tend to dissolve! We recommend that instructors test their equipment with the solvent that they intend to use before the start of the laboratory class.

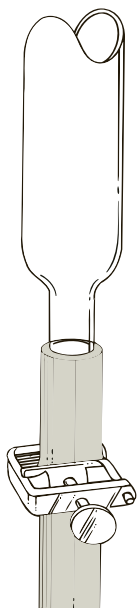


Figure 19.9
Tubing with screw clamp to regulate solvent flow on a chromatography column.

(Erlenmeyer flask) by adding the solid adsorbent, a little at a time, to a quantity of the solvent. This order of addition (adsorbent added to solvent) should be followed strictly because the adsorbent solvates and liberates heat. If the solvent is added to the adsorbent, it may boil away almost as fast as it is added due to heat evolved. This will be especially true if ether or another low-boiling solvent is used. When this happens, the final mixture will be uneven and lumpy. Enough adsorbent is added to the solvent, and mixed by swirling the container, to form a thick but flowing slurry. The container should be swirled until the mixture is homogeneous and relatively free of entrapped air bubbles.

For a standard-sized column, the procedure is as follows. When the slurry has been prepared, the column is filled about half full with solvent, and the stopcock is opened to allow solvent to drain slowly into a large beaker. The slurry is mixed by swirling and is then poured in portions into the top of the draining column (a wide-necked funnel may be useful here). Be sure to swirl the slurry thoroughly before each addition to the column. The column is tapped constantly and *gently* on the side during the pouring operation with the fingers or with a pencil fitted with a rubber stopper. A short piece of large-diameter pressure tubing may also be used for tapping. The tapping promotes even settling and mixing and gives an evenly packed column free of air bubbles. Tapping is continued until all the material has settled, showing a well-defined level at the top of the column. Solvent from the collecting beaker may be re-added to the slurry if it becomes too thick to be poured into the column at one time. In fact, the collected solvent should be cycled through the column several times to ensure that settling is complete and that the column is firmly packed. The downward flow of solvent tends to compact the adsorbent. You should take care never to let the column “run dry” during packing. There should always be solvent on top of the adsorbent column.

Dry Pack Method 1. In the first of the dry pack methods introduced here, the column is filled with solvent and allowed to drain *slowly*. The dry adsorbent is added, a little at a time, while the column is tapped gently with a pencil, finger, or glass rod.

Semimicroscale Columns. The procedure to fill a commercial semimicroscale column is essentially the same as that used to fill a Pasteur pipette (Section 19.6). The commercial column has the advantage that it is much easier to control the flow of solvent from the column during the filling process, because the stopcock can be adjusted appropriately. It is not necessary to use a cotton plug or to deposit a layer of sand before adding the adsorbent. The presence of the fritted disc at the base of the column acts to prevent adsorbent from escaping from the column.

Macroscale Columns. A plug of cotton is placed at the base of the column, and an even layer of sand is formed on top. The column is filled about half full with solvent, and the solid adsorbent is added carefully from a beaker while the solvent is allowed to flow slowly from the column. As the solid is added, the column is tapped as described for the slurry method in order to ensure that the column is packed evenly. When the column has the desired length, no more adsorbent is added. This method also produces an evenly packed column. Solvent should be cycled through this column (for macroscale applications) several times before each use. The same portion of solvent that has drained from the column during the packing is used to cycle through the column.

Dry Pack Method 2. In this method, the column is filled with dry adsorbent without any solvent. When the desired amount of adsorbent has been added, solvent is allowed to percolate through the column.

Semimicroscale Columns. The Dry Pack Method 2 is similar to that described for Pasteur pipettes (Section 19.6), except that the plug of cotton is not required. The flow rate of solvent through the column can be controlled using the stopcock, which is part of the column assembly (see Figure 19.8).

Macroscopic Columns. Macroscopic columns can also be packed by a dry pack method that is similar to the microscale methods described in Section 19.6. The disadvantages described for the microscale method also apply to the macroscopic method. This method is not recommended for use with silica gel or alumina, because the combination leads to uneven packing, air bubbles, and cracking, especially if a solvent that has a highly exothermic heat of solvation is used.

19.8 Applying the Sample to the Column

The solvent (or solvent mixture) used to pack the column is normally the least polar elution solvent that can be used during chromatography. The compounds to be chromatographed are not highly soluble in the solvent. If they were, they would probably have a greater affinity for the solvent than for the adsorbent and would pass right through the column without equilibrating with the stationary phase.

The first elution solvent, however, is generally not a good solvent to use in preparing the sample to be placed on the column. Because the compounds are not highly soluble in nonpolar solvents, it takes a large amount of the initial solvent to dissolve the compounds, and it is difficult to get the mixture to form a narrow band on top of the column. A narrow band is ideal for an optimum separation of components. For the best separation, therefore, the compound is applied to the top of the column undiluted if it is a liquid or in a *very small* amount of highly polar solvent if it is a solid. Water must not be used to dissolve the initial sample being chromatographed, because it reacts with the column packing.

In adding the sample to the column, use the following procedure. Lower the solvent level to the top of the adsorbent column by draining the solvent from the column. Add the sample (either a pure liquid or a solution) to form a small layer on top of the adsorbent. A Pasteur pipette is convenient for adding the sample to the column. Take care not to disturb the surface of the adsorbent. This is best accomplished by touching the pipette to the inside of the glass column and slowly draining it so as to allow the sample to spread into a thin film, which slowly descends to cover the entire adsorbent surface. Drain the pipette close to the surface of the adsorbent. When all the sample has been added, drain this small layer of liquid into the column until the top surface of the column *just begins* to dry. Then add a small layer of the chromatographic solvent carefully with a Pasteur pipette, again being careful not to disturb the surface. Drain this small layer of solvent into the column until the top surface of the column just dries. Add another small layer of fresh solvent, if necessary, and repeat the process until it is clear that the sample is strongly adsorbed on the top of the column. If the sample is colored and the fresh layer of solvent acquires some of this color, the sample has not been properly adsorbed. Once the sample has been properly applied, you can protect the level surface of the adsorbent by carefully filling the top of the column with solvent and sprinkling clean, white sand into the column so as to form a small protective layer on top of the adsorbent. For microscale applications, this layer of sand is not required.

Separations are often better if the sample is allowed to stand a short time on the column before elution. This allows a true equilibrium to be established. In columns that stand for too long, however, the adsorbent often compacts or even swells, and the flow can become annoyingly slow. Diffusion of the sample to widen the bands also becomes a problem if a column is allowed to stand over an extended period. For small-scale chromatography, using Pasteur pipettes, there is no stopcock, and it

is not possible to stop the flow. In this case, it is not considered necessary to allow the column to stand.

19.9 Elution Techniques

Solvents for analytical and preparative chromatography should be pure reagents. Commercial-grade solvents often contain small amounts of residue, which remains when the solvent is evaporated. For normal work, and for relatively easy separations that take only small amounts of solvent, the residue usually presents few problems. For large-scale work, commercial-grade solvents may have to be redistilled before use. This is especially true for hydrocarbon solvents, which tend to have more residue than other solvent types. Most of the experiments in this laboratory manual have been designed to avoid this particular problem.

One usually begins elution of the products with a nonpolar solvent, such as hexane or petroleum ether. The polarity of the elution solvent can be increased gradually by adding successively greater percentages of either ether or toluene (for instance, 1, 2, 5, 10, 15, 25, 50, 100%) or some other solvent of greater solvent power (polarity) than hexane. The transition from one solvent to another should not be too rapid in most solvent changes. If the two solvents to be changed differ greatly in their heats of solvation in binding to the adsorbent, enough heat can be generated to crack the column. Ether is especially troublesome in this respect because it has both a low boiling point and a relatively high heat of solvation. Most organic compounds can be separated on silica gel or alumina using hexane–ether or hexane–toluene combinations for elution and following these by pure methylene chloride. Solvents of greater polarity are usually avoided for the various reasons mentioned previously. In microscale work, the usual procedure is to use only one solvent for the chromatography.

The flow of solvent through the column should not be too rapid or the solutes will not have time to equilibrate with the adsorbent as they pass down the column. If the rate of flow is too low or stopped for a period, diffusion can become a problem—the solute band will diffuse, or spread out, in all directions. In either of these cases, separation will be poor. As a general rule (and only an approximate one), most macroscale columns are run with flow rates ranging from 5 to 50 drops of effluent per minute; a steady flow of solvent is usually avoided. Microscale columns made from Pasteur pipettes do not have a means of controlling the solvent flow rate, but commercial microscale columns are equipped with stopcocks. The solvent flow rate in this type of column can be adjusted in a manner similar to that used with larger columns. To avoid diffusion of the bands, do not stop the column and do not set it aside overnight.

19.10 Reservoirs

For microscale chromatography, the portion of the Pasteur pipette above the adsorbent is used as a reservoir of solvent. Fresh solvent, as needed, is added by means of another Pasteur pipette. When it is necessary to change solvent, the new solvent is also added in this manner. In some cases, the chromatography may proceed too slowly; the rate of solvent flow can be accelerated by attaching a rubber dropper bulb to the top of the Pasteur pipette column and squeezing *gently*. The additional air pressure forces the solvent through the column more rapidly. If this technique is used, however, care must be taken to remove the rubber bulb from the column before releasing it. Otherwise, air may be drawn up through the bottom of the column, destroying the column packing.

When large quantities of solvent are used in a chromatographic separation, it is often convenient to use a solvent reservoir to forestall having to add small portions of fresh solvent continually. The simplest type of reservoir, a feature of many columns, is created by fusing the top of the column to a round-bottom flask (Figure 19.10A). If

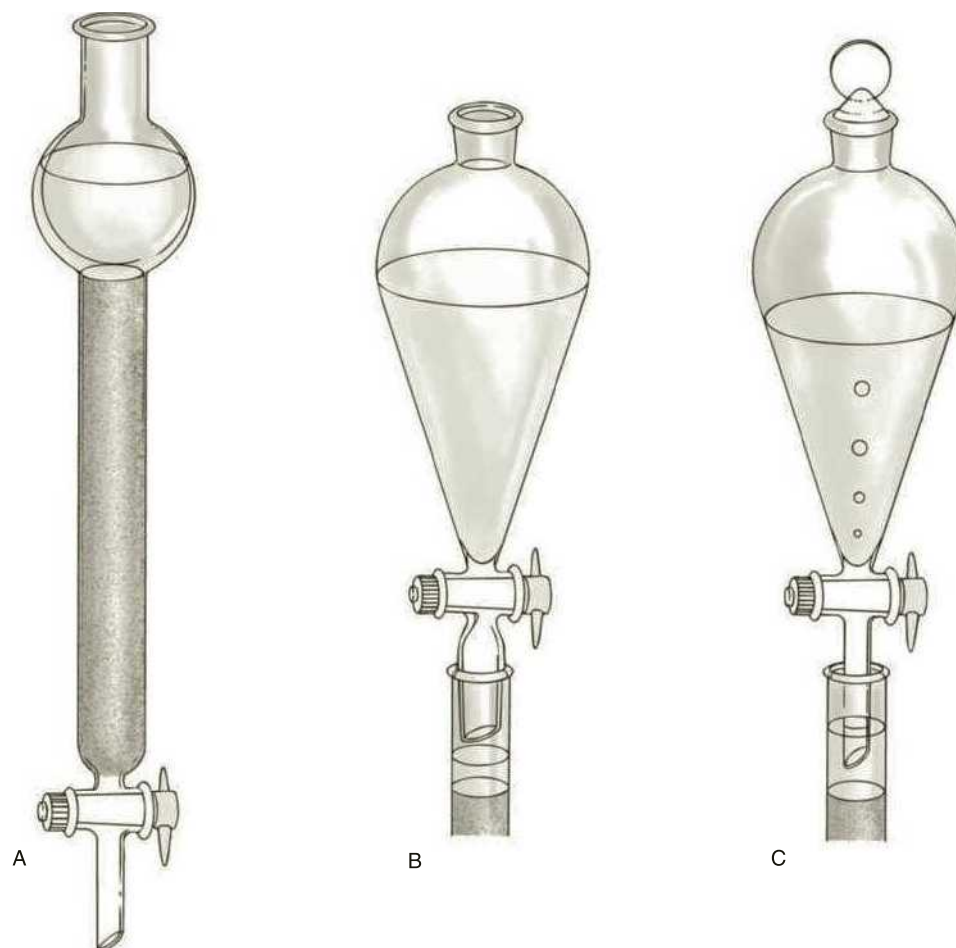


Figure 19.10
Various types of solvent-reservoir arrangements for chromatographic columns.

the column has a standard-taper joint at its top, a reservoir can be created by joining a standard-taper separatory funnel to the column (Figure 19.10B). In this arrangement, the stopcock is left open, and no stopper is placed in the top of the separatory funnel. A third common arrangement is shown in Figure 19.10C. A separatory funnel is filled with solvent; its stopper is wetted with solvent and put *firmly* in place. The funnel is inserted into the empty filling space at the top of the chromatographic column, and the stopcock is opened. Solvent flows out of the funnel, filling the space at the top of the column until the solvent level is well above the outlet of the separatory funnel. As solvent drains from the column, this arrangement automatically refills the space at the top of the column by allowing air to enter through the stem of the separatory funnel. Some microscale columns, such as that shown in Figure 19.8, are equipped with a solvent reservoir that fits onto the top of the column. It functions just like the reservoirs described in this section.

19.11 Monitoring the Column

It is a happy instance when the compounds to be separated are colored. The separation can then be followed visually and the various bands collected separately as they elute from the column. For the majority of organic compounds, however, this lucky circumstance does not exist, and other methods must be used to determine the positions of the bands. The most common method of following a separation of colorless compounds is to collect **fractions** of constant volume in preweighed

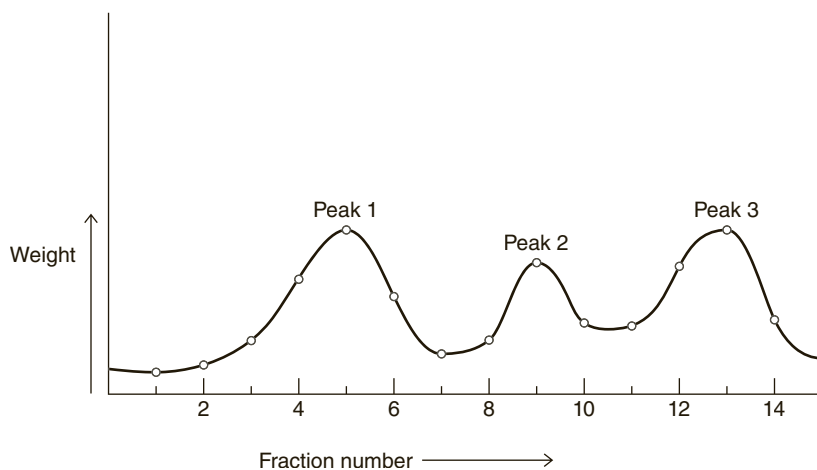


Figure 19.11
Typical elution graph.

flasks, to evaporate the solvent from each fraction, and to reweigh the flask plus any residue. A plot of fraction number versus the weight of the residues after evaporation of solvent gives a plot like that in Figure 19.11. Clearly, fractions 2 through 7 (Peak 1) may be combined as a single compound, and so can fractions 8 through 11 (Peak 2) and 12 through 15 (Peak 3). The size of the fractions collected (1, 10, 100, or 500 mL) depends on the size of the column and the ease of separation.

Another common method of monitoring the column is to mix an inorganic phosphor into the adsorbent used to pack the column. When the column is illuminated with an ultraviolet light, the adsorbent treated in this way fluoresces. However, many solutes have the ability to **quench** the fluorescence of the indicator phosphor. In areas in which solutes are present, the adsorbent does not fluoresce, and a dark band is visible. In this type of column, the separation can also be followed visually.

Thin-layer chromatography is often used to monitor a column. This method is described in Technique 20 (Section 20.10). Several sophisticated instrumental and spectroscopic methods, which we shall not detail, can also monitor a chromatographic separation.

19.12 Tailing

When a single solvent is used for elution, an elution curve (weight versus fraction) like that shown as a solid line in Figure 19.12 is often observed. An ideal elution curve is shown by dashed lines. In the nonideal curve, the compound is said to be **tailing**. Tailing can interfere with the beginning of a curve or a peak of a second component and lead to a poor separation. One way to avoid this is to increase the polarity of the solvent constantly while eluting. In this way, at the tail of the peak, where the solvent polarity is increasing, the compound will move slightly faster than at the front and allow the tail to squeeze forward, forming a more nearly ideal band.

19.13 Recovering the Separated Compounds

In recovering each of the separated compounds of a chromatographic separation when they are solids, the various correct fractions are combined and evaporated. If the combined fractions contain sufficient material, they may be purified by recrystallization. If the compounds are liquids, the correct fractions are combined, and the solvent is evaporated. If sufficient material has been collected, liquid samples can be purified by distillation. The combination of chromatography–crystallization or

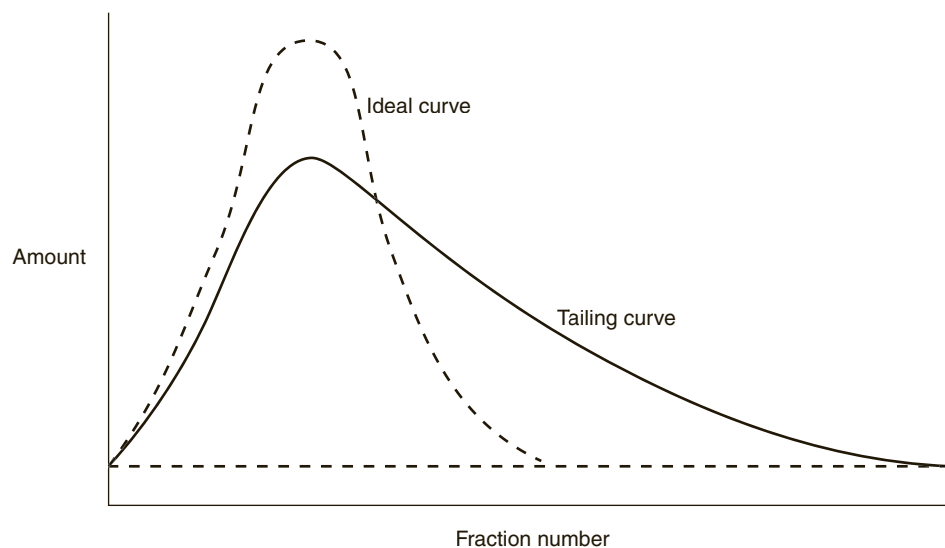


Figure 19.12
Elution curves: One ideal and one that “tails.”

chromatography–distillation usually yields very pure compounds. For microscale applications, the amount of sample collected is too small to allow a purification by crystallization or distillation. The samples that are obtained after the solvent has been evaporated are considered to be sufficiently pure, and no additional purification is attempted.

19.14 Decolorization by Column Chromatography

A common outcome of organic reactions is the formation of a product that is contaminated by highly colored impurities. Very often these impurities are highly polar, and they have a high molecular weight, as well as being colored. The purification of the desired product requires that these impurities be removed. Section 11.7 of Technique 11 details methods of decolorizing an organic product. In most cases, these methods involve the use of a form of activated charcoal, or Norit.

An alternative, which is applied conveniently in microscale experiments, is to remove the colored impurity by column chromatography. Because of the polarity of the impurities, the colored components are strongly adsorbed on the stationary phase of the column, and the less polar desired product passes through the column and is collected.

Microscale decolorization of a solution on a chromatography column requires that a column be prepared in a Pasteur pipette, using either alumina or silica gel as the adsorbent (Section 19.6). The sample to be decolorized is diluted to the point where crystallization within the column will not take place, and it is then passed through the column in the usual manner. The desired compound is collected as it exits the column, and the excess solvent is removed by evaporation (Technique 7, Section 7.10).

19.15 Gel Chromatography

The stationary phase in gel chromatography consists of a cross-linked polymeric material. Molecules are separated according to their *size* by their ability to penetrate a sieve-like structure. Molecules permeate the porous stationary phase as they move down the column. Small molecules penetrate the porous structure more easily than large ones. Thus, the large molecules move through the column faster than the smaller ones and elute first. The separation of molecules by gel chromatography is depicted in Figure 19.13. With adsorption chromatography using materials

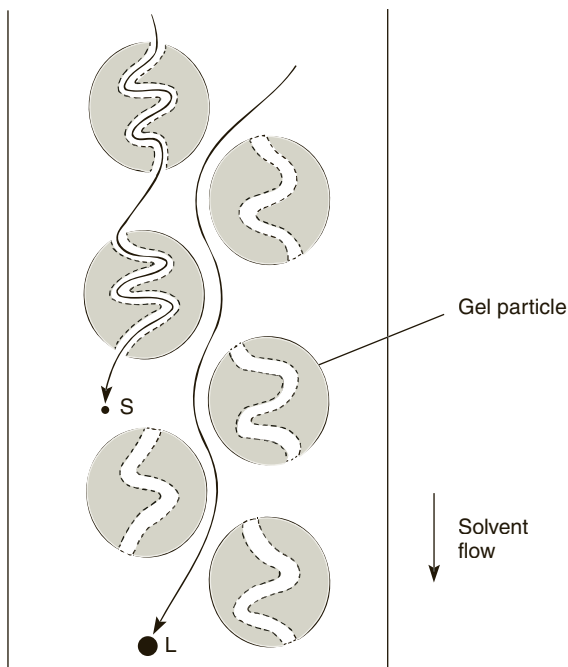


Figure 19.13
Gel chromatography: Comparison of the paths of large (L) and small (S) molecules through the column during the same interval.

such as alumina or silica, the order is usually the reverse. Small molecules (of low molecular weight) pass through the column *faster* than large molecules (of high molecular weight) because large molecules are more strongly attracted to the polar stationary phase.


Equivalent terms used by chemists for the gel-chromatography technique are **gel filtration** (biochemistry term), **gel-permeation chromatography** (polymer chemistry term), and **molecular sieve chromatography**. **Size-exclusion chromatography** is a general term for the technique, and it is perhaps the most descriptive term for what occurs on a molecular level.

Sephadex is one of the most popular materials for gel chromatography. It is widely used by biochemists for separating proteins, nucleic acids, enzymes, and carbohydrates. Most often, water or aqueous solutions of buffers are used as the moving phase. Chemically, Sephadex is a polymeric carbohydrate that has been cross-linked. The degree of cross-linking determines the size of the “holes” in the polymer matrix. In addition, the hydroxyl groups on the polymer can adsorb water, which causes the material to swell. As it expands, “holes” are created in the matrix. Several different gels are available from manufacturers, each with its own set of characteristics.

For example, a typical Sephadex gel, such as G-75, can separate molecules in the molecular weight (MW) range 3000 to 70,000. Assume for the moment that one has a four-component mixture containing compounds with molecular weights of 10,000, 20,000, 50,000, and 100,000. The 100,000-MW compound would pass through the column first because it cannot penetrate the polymer matrix. The 50,000-, 20,000-, and 10,000-MW compounds penetrate the matrix to various degrees and would be separated. The molecules would elute in the order given (decreasing order of molecular weights). The gel separates on the basis of molecular size and configuration, rather than molecular weight.

Sephadex LH-20 has been developed for nonaqueous solvents. Some of the hydroxyl groups have been alkylated, and thus the material can swell under both aqueous and nonaqueous conditions (it now has “organic” character). This material can be used with several organic solvents, such as alcohol, acetone, methylene chloride, and aromatic hydrocarbons.

Another type of gel is based on a polyacrylamide structure (Bio-Gel P and Poly-Sep AA). A portion of a polyacrylamide chain is shown here:

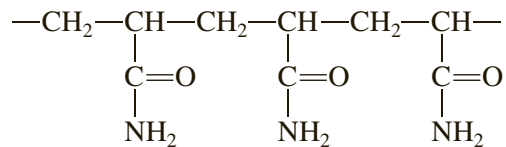

Gels of this type can also be used in water and some polar organic solvents. They tend to be more stable than Sephadex, especially under acidic conditions. Polyacrylamides can be used for many biochemical applications involving macromolecules. For separating synthetic polymers, cross-linked polystyrene beads (copolymer of styrene and divinylbenzene) find common application. Again, the beads are swollen before use. Common organic solvents can be used to elute the

polymers. As with other gels, the higher-molecular-weight compounds elute before the lower-molecular-weight compounds.

19.16 Flash Chromatography

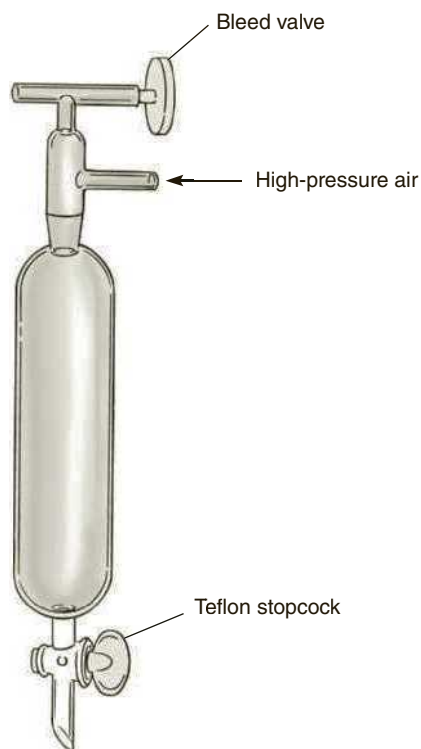


Figure 19.14
Apparatus for flash chromatography.

One of the drawbacks to column chromatography is that, for large-scale preparative separations, the time required to complete a separation may be very long.

Furthermore, the resolution that is possible for a particular experiment tends to deteriorate as the time for the experiment grows longer. This latter effect arises because the bands of compounds that move very slowly through a column tend to “tail.”

A technique that can be useful in overcoming these problems has been developed. This technique, called **flash chromatography**, is actually a simple modification of ordinary column chromatography. In flash chromatography, the adsorbent is packed into a relatively short glass column, and air pressure is used to force the solvent through the adsorbent.

The apparatus used for flash chromatography is shown in Figure 19.14. The glass column is fitted with a Teflon stopcock at the bottom to control the flow rate of solvent. A plug of glass wool is placed in the bottom of the column to act as a support for the adsorbent. A layer of sand may also be added on top of the glass wool. The column is filled with adsorbent using the dry pack method. When the column has been filled, a fitting is attached to the top of the column, and the entire apparatus is connected to a source of high-pressure air or nitrogen. The fitting is designed so that the pressure applied to the top of the column can be adjusted precisely. The source of the high-pressure air is often a specially adapted air pump.

A typical column would use silica gel adsorbent (particle size = 40 to 63 μm) packed to a height of 5 inches in a glass column of 20-mm diameter. The pressure applied to the column would be adjusted to achieve a solvent flow rate such that the solvent level in the column would decrease by about 2 in./min. This system would be appropriate to separate the components of a 250-mg sample.

The high-pressure air forces the solvent through the column of adsorbent at a rate that is much greater than what would be achieved if the solvent flowed through the column under the force of gravity. Because the solvent is caused to flow faster, the time required for substances to pass through the column is reduced. By itself, simply applying air pressure to the column might reduce the clarity of the separation because the components of the mixture would not have time to establish themselves into distinctly separate bands. However, in flash chromatography, you can use a much finer adsorbent than would be used in ordinary chromatography. With a much smaller particle size for the adsorbent, the surface area is increased, and the resolution possible thereby improves.

A simple variation on this idea does not use air pressure. Instead, the lower end of the column is inserted into a stopper, which is fitted into the top of a suction flask. Vacuum is applied to the system, and the vacuum acts to draw the solvent through the adsorbent column. The overall effect of this variation is similar to that obtained when air pressure is applied to the top of the column.

REFERENCES

- Deyl, Z.; Macek, K.; Janák, J. *Liquid Column Chromatography*. Elsevier: Amsterdam, 1975.
- Heftmann, E. *Chromatography*, 3rd ed. Van Nostrand Reinhold: New York, 1975.
- Jacobson, B. M. An Inexpensive Way to Do Flash Chromatography. *Journal of Chemical Education*, 65 (May 1988): 459.
- Still, W. C.; Kahn, M.; Mitra, A Rapid Chromatographic Technique for Preparative Separations with Moderate Resolution. *Journal of Organic Chemistry*, 43 (1978): 2923.

PROBLEMS

1. A sample was placed on a chromatography column. Methylene chloride was used as the eluting solvent. No separation of the components in the sample was observed. What must have been happening during this experiment? How would you change the experiment to overcome this problem?
2. You are about to purify an impure sample of naphthalene by column chromatography. What solvent should you use to elute the sample?
3. Consider a sample that is a mixture composed of biphenyl, benzoic acid, and benzyl alcohol. Predict the order of elution of the components in this mixture. Assume that the chromatography uses a silica column, and the solvent system is based on cyclohexane, with an increasing proportion of methylene chloride being added as a function of time.
4. An orange compound was added to the top of a chromatography column. Solvent was added immediately, with the result that the entire volume of solvent in the solvent reservoir turned orange. No separation could be obtained from the chromatography experiment. What went wrong?
5. A yellow compound, dissolved in methylene chloride, is added to a chromatography column. The elution is begun using petroleum ether as the solvent. After 6 L of solvent have passed through the column, the yellow band still has not traveled down the column appreciably. What should be done to make this experiment work better?
6. You have 0.50 g of a mixture that you wish to purify by column chromatography. How much adsorbent should you use to pack the column? Estimate the appropriate column diameter and height.
7. In a particular sample, you wish to collect the component with the *highest* molecular weight as the *first* fraction. What chromatographic technique should you use?
8. A colored band shows an excessive amount of tailing as it passes through the column. What can you do to rectify this problem?
9. How would you monitor the progress of a column chromatography when the sample is colorless? Describe at least two methods.

Thin-Layer Chromatography

Thin-layer chromatography (TLC) is a very important technique for the rapid separation and qualitative analysis of small amounts of material. It is ideally suited for the analysis of mixtures and reaction products in microscale experiments. The technique is closely related to column chromatography. In fact, TLC can simply be considered column chromatography in reverse, with the solvent ascending the adsorbent, rather than descending. Because of this close relationship to column chromatography, and because the principles governing the two techniques are similar, Technique 19, on column chromatography, should be read first.

20.1 Principles of Thin-Layer Chromatography

Like column chromatography, TLC is a solid–liquid partitioning technique. However, the moving liquid phase is not allowed to percolate down the adsorbent; it is caused to ascend a thin layer of adsorbent coated onto a backing support. The most typical backing is a glass plate, but other materials are also used. A thin layer of the adsorbent is spread onto the plate and allowed to dry. A coated and dried plate of glass is called a **thin-layer plate** or a **thin-layer slide**. (The reference to *slide* comes about because microscope slides are often used to prepare small thin-layer plates.) When a thin-layer plate is placed upright in a vessel that contains a shallow layer of solvent, the solvent ascends the layer of adsorbent on the plate by capillary action.

In TLC, the sample is applied to the plate before the solvent is allowed to ascend the adsorbent layer. The sample is usually applied as a small spot near the base of the plate; this technique is often referred to as **spotting**. The plate is spotted by repeated applications of a sample solution from a small capillary pipette. When the filled pipette touches the plate, capillary action delivers its contents to the plate, and a small spot is formed.

As the solvent ascends the plate, the sample is partitioned between the moving liquid phase and the stationary solid phase. During this process, you are **developing**, or **running**, the thin-layer plate. In development, the various components in the applied mixture are separated. The separation is based on the many equilibrations the solutes experience between the moving and the stationary phases. (The nature of these equilibrations was thoroughly discussed in Technique 19, Sections 19.2 and 19.3.) As in column chromatography, the least polar substances advance faster than the most polar substances. A separation results from the differences in the rates at which the individual components of the mixture advance upward on the plate. When many substances are present in a mixture, each has its own characteristic solubility and adsorptivity properties, depending on the functional groups in its structure. In general, the stationary phase is strongly polar and strongly binds polar substances. The moving liquid phase is usually less polar than the adsorbent and most easily dissolves substances that are less polar or even nonpolar. Thus, substances that are the most polar travel slowly upward, or not at all, and nonpolar substances travel more rapidly if the solvent is sufficiently nonpolar.

When the thin-layer plate has been developed, it is removed from the developing tank and allowed to dry until it is free of solvent. If the mixture that was originally spotted on the plate was separated, there will be a vertical series of spots on the plate. Each spot corresponds to a separate component or compound from the original mixture. If the components of the mixture are colored substances, the various spots will be clearly visible after development. More often, however, the “spots” will not be visible because they correspond to colorless substances. If spots are not apparent, they can be made visible only if a **visualization method** is used. Often, spots can be seen when the thin-layer plate is held under ultraviolet light; the ultraviolet lamp is a common visualization method. Also common is the use of iodine vapor. The plates are placed in a chamber containing iodine crystals and left to stand for a short time. The iodine reacts with the various compounds adsorbed on the plate to give colored complexes that are clearly visible. Because iodine often changes the compounds by reaction, the components of the mixture cannot be recovered from the plate when the iodine method is used. (Other methods of visualization are discussed in Section 20.7.)

20.2 Commercially Prepared TLC Plates

The most convenient type of TLC plate is prepared commercially and sold in a ready-to-use form. Many manufacturers supply glass plates precoated with a durable layer of silica gel or alumina. More conveniently, plates are also available that have either a flexible plastic backing or an aluminum backing. The most common types of commercial TLC plates are composed of plastic sheets that are coated with silica gel and polyacrylic acid, which serves as a binder. A fluorescent indicator may be mixed with the silica gel. The indicator renders the spots due to the presence of compounds in the sample visible under ultraviolet light (see Section 20.7). Although these plates are relatively expensive when compared with plates prepared in the laboratory, they are far more convenient to use, and they provide more consistent results. The plates are manufactured quite uniformly. Because the plastic backing is flexible, they have the additional advantage that the coating does not flake off the plates easily. The plastic sheets (usually 8 in. by 8 in. square) can also be cut with a pair of scissors or paper cutter to whatever size may be required.

If the package of commercially prepared TLC plates has been opened previously, or if the plates have not been purchased recently, they should be dried before use. Dry the plates by placing them in an oven at 100°C for 30 minutes, and store them in a desiccator until they are to be used.

20.3 Preparation of Thin-Layer Slides and Plates

Commercially prepared plates (Section 20.2) are the most convenient to use, but if you must prepare your own slides or plates, this section gives the directions for preparing them. The two adsorbent materials used most often for TLC are alumina G (aluminum oxide) and silica gel G (silicic acid). The G designation stands for gypsum (calcium sulfate). Calcined gypsum, $\text{CaSO}_4 \cdot \frac{1}{2}\text{H}_2\text{O}$, is better known as plaster of Paris. When exposed to water or moisture, gypsum sets in a rigid mass, $\text{CaSO}_4 \cdot 2\text{H}_2\text{O}$, which binds the adsorbent together and to the glass plates used as a backing support. In the adsorbents used for TLC, about 10–13% by weight of gypsum is added as a binder. The adsorbent materials are otherwise like those used in column chromatography; the adsorbents used in column chromatography have a larger particle size, however. The material for thin-layer work is a fine powder. The small particle size, along with the added gypsum, makes it impossible to use silica gel G or alumina G for column work. In a

column, these adsorbents generally set so rigidly that solvent virtually stops flowing through the column.

A. Microscope Slide TLC Plates

For qualitative work such as identifying the number of components in a mixture or trying to establish that two compounds are identical, small TLC plates made from microscope slides are especially convenient. Coated microscope slides are easily made by dipping the slides into a container holding a slurry of the adsorbent material. Although numerous solvents can be used to prepare a slurry, methylene chloride is probably the best choice. It has the two advantages of low boiling point (40°C) and inability to cause the adsorbent to set or form lumps. The low boiling point means that it is not necessary to dry the coated slides in an oven. Its inability to cause the gypsum binder to set means that slurries made with it are stable for several days. The layer of adsorbent formed is fragile, however, and must be treated carefully. For this reason, some people prefer to add a small amount of methanol to the methylene chloride to enable the gypsum to set more firmly. The methanol solvates the calcium sulfate much as water does. More durable plates can be made by dipping plates into a slurry prepared from water. These plates must be oven-dried before use. Also, a slurry prepared from water must be used soon after its preparation. If it is not, it begins to set and form lumps. Thus, an aqueous slurry must be prepared immediately before use; it cannot be used after it has stood for any length of time. For microscope slides, a slurry of silica gel G in methylene chloride is not only convenient but also adequate for most purposes.

Preparing the Slurry. The slurry is most conveniently prepared in a 4-oz wide-mouthed screw-cap jar. About 3 mL of methylene chloride are required for each gram of silica gel G. For a smooth slurry without lumps, the silica gel should be added to the solvent while the mixture is being either stirred or swirled. Adding solvent to the adsorbent usually causes lumps to form in the mixture. When the addition is complete, the cap should be placed on the jar tightly and the jar shaken vigorously to ensure thorough mixing. The slurry may be stored in the tightly capped jar until it is to be used. More methylene chloride may have to be added to replace evaporation losses.

CAUTION



Avoid breathing silica dust or methylene chloride, prepare and use the slurry in a hood, and avoid getting methylene chloride or the slurry mixture on your skin. Wear protective gloves.

Preparing the Slides. If new microscope slides are available, you can use them without any special treatment. However, it is more economical to reuse or recycle microscope slides. Wash the slides with soap and water, rinse with water, and then rinse with 50% aqueous methanol. Allow the plates to dry thoroughly on paper towels. They should be handled by the edges because fingerprints on the plate surface will make it difficult for the adsorbent to bind to the glass.

Coating the Slides. The slides are coated with adsorbent by dipping them into the container of slurry. You can coat two slides simultaneously by sandwiching them together before dipping them in the slurry.

NOTE: Perform the coating operation under a hood.

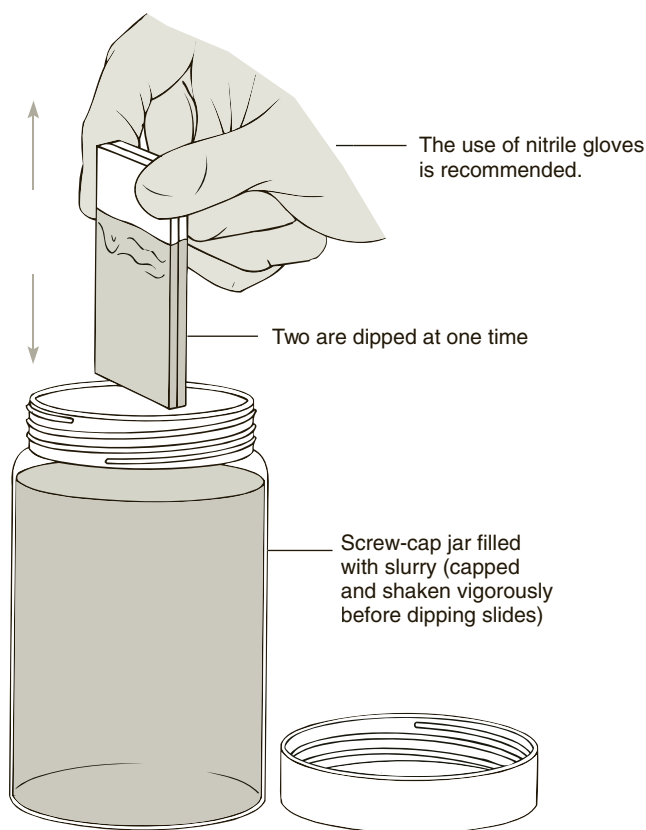


Figure 20.1
Dipping slides to coat them.

Shake the slurry vigorously just before dipping the slides. Because the slurry settles on standing, it should be mixed in this way before each set of slides is dipped. The depth of the slurry in the jar should be about 3 inches, and the plates should be dipped into the slurry until only about 0.25 inches at the top remains uncoated. The dipping operation should be performed smoothly. The plates may be held at the top (see Figure 20.1), where they will not be coated. They are dipped into the slurry and withdrawn with a slow and steady motion. The dipping operation takes about 2 seconds. Some practice may be required to get the correct timing. After dipping, replace the cap on the jar and hold the plates for a minute until most of the solvent has evaporated. Separate the plates and place them on paper towels to complete the drying. The plates should have an even coating; there should be no streaks and no thin spots where glass shows through the adsorbent. The plates should not have a thick and lumpy coating.

Two conditions cause thin and streaked plates. First, the slurry may not have been mixed thoroughly before the dipping operation; the adsorbent might then have settled to the bottom of the jar, and the thin slurry at the top would not have coated the slides properly. Second, the slurry simply may

not have been thick enough; more silica gel G must then be added to the slurry until the consistency is correct. If the slurry is too thick, the coating on the plates will be thick, uneven, and lumpy. To correct this, dilute the slurry with enough solvent to achieve the proper consistency.

Plates with an unsatisfactory coating may be wiped clean with a paper towel and redipped. Take care to handle the plates only from the top or by the sides to avoid getting fingerprints on the glass surface.

B. Larger Thin-Layer Plates

For separations involving large amounts of material, or for difficult separations, it may be necessary to use larger thin-layer plates. Plates with dimensions up to 20–25 cm² are common. With larger plates, it is desirable to have a more durable coating, and a water slurry of the adsorbent should be used to prepare them. If silica gel is used, the slurry should be prepared in the ratio of about 1 g silica gel G to each 2 mL of water. The glass plate used for the thin-layer plate should be washed, dried, and placed on a sheet of newspaper. Place two strips of masking tape along two edges of the plate. Use more than one layer of masking tape if a thicker coating is desired on the plate. A slurry is prepared, shaken well, and poured along one of the untaped edges of the plate.

NOTE: Observe the same precautions stated for Microscope Slide TLC Plates.

A heavy piece of glass rod, long enough to span the taped edges, is used to level and spread the slurry over the plate. While the rod is resting on the tape, it is pushed

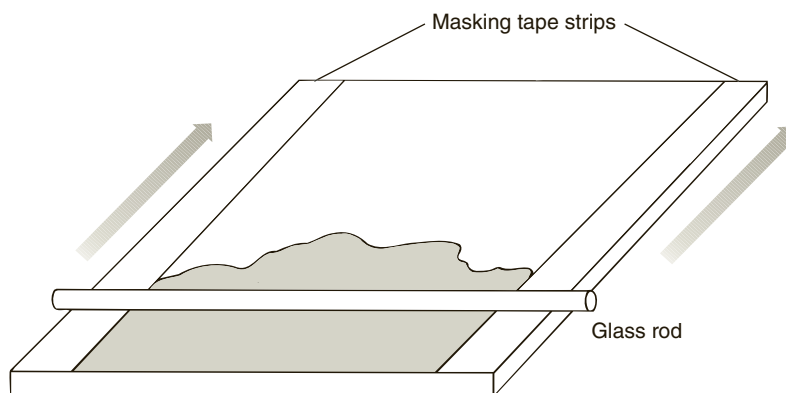


Figure 20.2
Preparing a large plate.

along the plate from the end at which the slurry was poured toward the opposite end of the plate. This is illustrated in Figure 20.2. After the slurry is spread, the masking tape strips are removed, and the plates are dried in a 110°C oven for about 1 hour. Plates of 10–25 cm² are easily prepared by this method. Larger plates present more difficulties. Many laboratories have a commercially manufactured spreading machine that makes the entire operation simpler.

20.4 Sample Application: Spotting the Plates

Preparing a Micropipette

To apply the sample that is to be separated to the thin-layer plate, one uses a micropipette. A micropipette is easily made from a short length of thin-walled capillary tubing like that used for melting-point determinations, but open at both ends. The capillary tubing is heated at its midpoint with a microburner and rotated until it is soft. When the tubing is soft, the heated portion of the tubing is drawn out until a constricted portion of tubing 4–5 cm long is formed. After cooling, the constricted portion of tubing is scored at its center with a file or scorer and broken. The two halves yield two capillary micropipettes. Figure 20.3 shows how to make such pipettes.

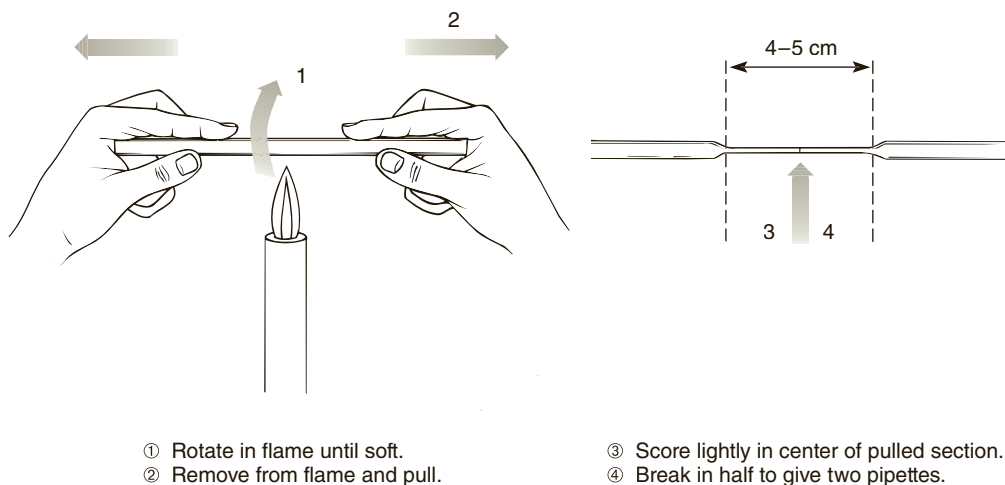


Figure 20.3
Construction of two capillary micropipettes.

Spotting the Plate

To apply a sample to the plate, begin by placing about 1 mg of a solid test substance, or one drop of a liquid test substance, in a small container such as a watch glass or a test

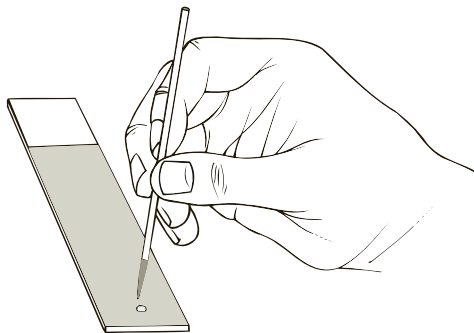


Figure 20.4
Spotting the plate with a drawn capillary pipette.

tube. Dissolve the sample in a few drops of a volatile solvent. Acetone or methylene chloride is usually a suitable solvent. If a solution is to be tested, it can often be used directly. The small capillary pipette, prepared as described, is filled by dipping the pulled end into the solution to be examined. Capillary action fills the pipette. Empty the pipette by touching it *lightly* to the thin-layer plate at a point about 1 cm from the bottom (Figure 20.4). The spot must be high enough that it does not dissolve in the developing solvent. It is important to touch the plate very lightly and not to gouge a hole in the adsorbent. When the pipette touches the plate, the solution is transferred to the plate as a small spot. The pipette should be touched to the plate *very briefly* and then removed. If the pipette is held to the plate, its entire contents will be delivered to the plate.

Only a small amount of material is needed. It is often helpful to blow gently on the plate as the sample is applied. This helps to keep the spot small by evaporating the solvent before it can spread out on the plate. The smaller the spot formed, the better the separation obtainable. If needed, additional material can be applied to the plate by repeating the spotting procedure. You should repeat the procedure with several small amounts, rather than apply one large amount. The solvent should be allowed to evaporate between applications. If the spot is not small (about 2 mm in diameter), a new plate should be prepared. The capillary pipette may be used several times if it is rinsed between uses. It is repeatedly dipped into a small portion of solvent to rinse it and touched to a paper towel to empty it.

As many as three spots may be applied to a microscope-slide TLC plate. Each spot should be about 1 cm from the bottom of the plate, and all spots should be evenly spaced, about 1 cm apart. The spots should be positioned at least 1 cm from the edges of the plate. Due to diffusion, spots often increase in diameter as the plate is developed. To keep spots containing different materials from merging and to avoid confusing the samples, do not place more than three spots on a single plate. Larger plates can accommodate many more samples.

20.5 Developing (Running) TLC Plates

Preparing a Development Chamber

A convenient development chamber for microscope-slide TLC plates can be made from a 4-oz wide-mouthed jar. An alternative development chamber can be constructed from a beaker, using aluminum foil to cover the opening. The inside of the jar or beaker should be lined with a piece of filter paper, cut so that it does not quite extend around the inside of the jar. A small vertical opening (2–3 cm) should be left in the filter paper for observing the development. Before development, the filter paper inside the jar or beaker should be thoroughly moistened with the development solvent. The solvent-saturated liner helps to keep the chamber saturated with solvent vapors, thereby speeding the development. Once the liner is saturated, the level of solvent in the bottom of the development chamber is adjusted to a depth of about 5 mm, and the chamber is capped (or covered with aluminum foil) and set aside until it is to be used. A correctly prepared development chamber (with slide in place) is shown in Figure 20.5.

Developing the TLC Plate

Once the spot has been applied to the thin-layer plate and the solvent has been selected (see Section 20.5), the plate is placed in the chamber for development. The

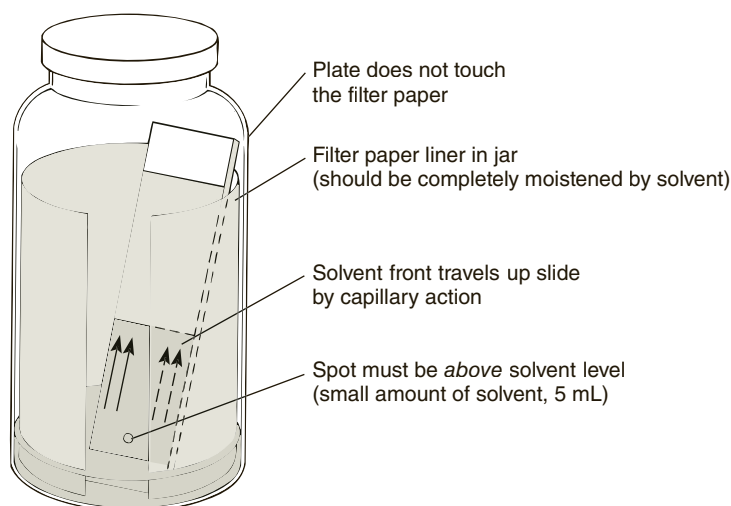


Figure 20.5
Development chamber with thin-layer plate undergoing development.

plate must be placed in the chamber carefully so that none of the coated portion touches the filter paper liner. In addition, the solvent level in the bottom of the chamber must not be above the spot that was applied to the plate, or the spotted material will dissolve in the pool of solvent instead of undergoing chromatography. Once the plate has been placed correctly, replace the cap on the developing chamber and wait for the solvent to advance up the plate by capillary action. This generally occurs rapidly, and you should watch carefully. As the solvent rises, the plate becomes visibly moist. When the solvent has advanced to within 5 mm of the end of the coated surface, the plate should be removed, and the position of the solvent front should be marked *immediately* by scoring the plate along the solvent line with a *pencil*. The solvent front must not be allowed to travel beyond the end of the coated surface. The plate should be removed before this happens. The solvent will not actually advance beyond the end of the plate, but spots allowed to stand on a completely moistened plate on which the solvent is not in motion expand by diffusion. Once the plate has dried, any visible spots should be outlined on the plate with a pencil. If no spots are apparent, a visualization method (Section 20.7) may be needed.

20.6 Choosing a Solvent for Development

The development solvent used depends on the materials to be separated. You may have to try several solvents before a satisfactory separation is achieved. Because microscope slides can be prepared and developed rapidly, an empirical choice is usually not hard to make. A solvent that causes all the spotted material to move with the solvent front is too polar. One that does not cause any of the material in the spot to move is not polar enough. As a guide to the relative polarity of solvents, consult Table 19.2 in Technique 19. Figure 20.6 shows three TLC plates run with different solvents. As can be seen, if the solvent is too polar, the spots tend to run near the top of the plate, and the separation is poor. If the solvent is not sufficiently polar, the spots do not travel very far; again, the separation is poor. The ideal solvent choice allows the spots to travel well up the plate but allows for a clean separation.

Methylene chloride and toluene are solvents of intermediate polarity and good choices for a wide variety of functional groups to be separated. For hydrocarbon materials, good first choices are hexane, petroleum ether (ligroin), or toluene. Hexane or petroleum ether with varying proportions of toluene or ether give

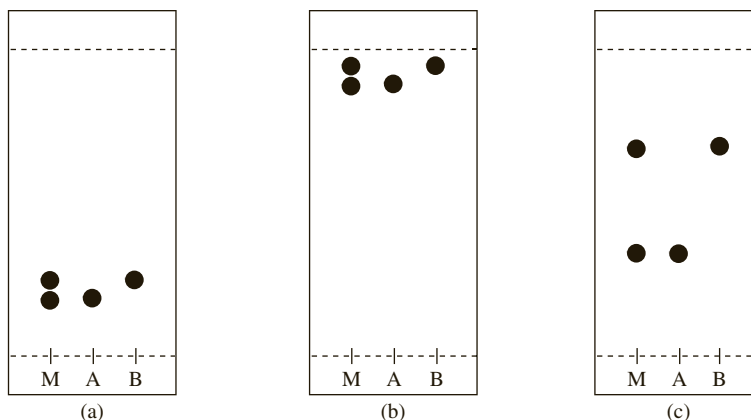


Figure 20.6

TLC plates showing the effects of solvent polarity on a separation.

The three spots on each plate are indicated by M for a two-component mixture, A for a standard sample of substance A, and B for a standard sample of substance B. (a) Solvent of low polarity: poor separation. (b) Solvent of high polarity: poor separation. (c) Solvent of intermediate polarity: clean separation.

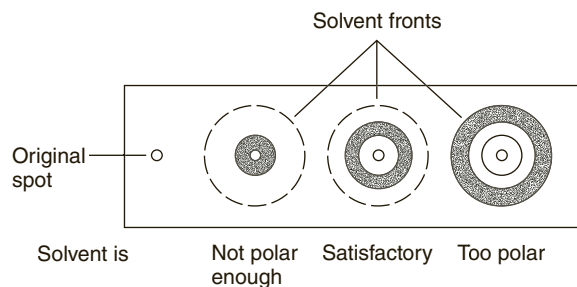


Figure 20.7

Concentric ring method of testing solvents.

solvent mixtures of moderate polarity that are useful for many common functional groups. Polar materials may require ethyl acetate, acetone, or methanol.

A rapid way to determine a good solvent is to apply several sample spots to a single plate. The spots should be placed a minimum of 1 cm apart. A capillary pipette is filled with a solvent and gently touched to one of the spots. The solvent expands outward in a circle. The solvent front should be marked with a pencil. A different solvent is applied to each spot. As the solvents expand outward, the spots expand as concentric rings. From the appearance of the rings, you can judge approximately the suitability of the solvent. Several types of behavior experienced with this method of testing are shown in Figure 20.7.

20.7 Visualization Methods

It is fortunate when the compounds separated by TLC are colored because the separation can be followed visually. More often than not, however, the compounds are colorless. In that case, the separated materials must be made visible by some reagent or some method that makes the separated compounds visible. Reagents that give rise to colored spots are called **visualization reagents**. Methods of viewing that make the spots apparent are **visualization methods**.

The visualization reagent used most often is iodine. Iodine reacts with many organic materials to form complexes that are either brown or yellow. In this visualization method, the developed and dried TLC plate is placed in a 4-oz wide-mouth

screw-cap jar along with a few crystals of iodine. The jar is capped and gently warmed on a steam bath or a hot plate at low heat. The jar fills with iodine vapors, and the spots begin to appear. When the spots are sufficiently intense, the plate is removed from the jar and the spots are outlined with a pencil. The spots are not permanent. Their appearance results from the formation of complexes the iodine makes with the organic substances. As the iodine sublimes off the plate, the spots fade. Hence, they should be marked immediately. Nearly all compounds except saturated hydrocarbons and alkyl halides form complexes with iodine. The intensities of the spots do not accurately indicate the amount of material present, except in the crudest way.

A more convenient variation on this method is to combine 50 grams of silica gel with about 200 mg of iodine crystals in an 8-oz. jar.¹ Rotate the jar periodically to disperse the iodine. Making certain that the TLC plate is completely dry, immerse it in the silica gel-iodine mixture and close the lid tightly. Rotate the jar to ensure that the TLC plate is completely covered. After one or two minutes, remove the plate from the jar using forceps. Be careful not to inhale the powder!

The second most common method of visualization is by an ultraviolet (UV) lamp. Under UV light, compounds often look like bright spots on the plate. This often suggests the structure of the compound: Certain types of compounds shine very brightly under UV light because they fluoresce.

Another method that provides good results involves adding a fluorescent indicator to the adsorbent used to coat the plates. A mixture of zinc and cadmium sulfides is often used. When treated in this way and held under UV light, the entire plate fluoresces. However, dark spots appear on the plate where the separated compounds are seen to quench this fluorescence.

In addition to the preceding methods, several chemical methods are available that either destroy or permanently alter the separated compounds through reaction. Many of these methods are specific for particular functional groups.

Alkyl halides can be visualized if a dilute solution of silver nitrate is sprayed on the plates. Silver halides are formed. These halides decompose if exposed to light, giving rise to dark spots (free silver) on the TLC plate.

Most organic functional groups can be made visible if they are charred with sulfuric acid. Concentrated sulfuric acid is sprayed on the plate, which is then heated in an oven at 110°C to complete the charring. Permanent spots are thus created.

Colored compounds can be prepared from colorless compounds by making derivatives before spotting them on the plate. An example of this is the preparation of 2,4-dinitrophenylhydrazones from aldehydes and ketones to produce yellow and orange compounds. You may also spray the 2,4-dinitrophenylhydrazine reagent on the plate after the ketones or aldehydes have separated. Red and yellow spots form where the compounds are located. Other examples of this method are using ferric chloride for visualizing phenols and using bromocresol green for detecting carboxylic acids. Chromium trioxide, potassium dichromate, and potassium permanganate can be used for visualizing compounds that are easily oxidized. *p*-Dimethylaminobenzaldehyde easily detects amines. Ninhydrin reacts with amino acids to make them visible. Numerous other methods and reagents available from various supply outlets are specific for certain types of functional groups. These visualize only the class of compounds of interest.

¹The authors are grateful to Prof. Amanda Murphy for suggesting this method. It is best to prepare the jar at least one day in advance of class period to allow the iodine to disperse evenly throughout the silica gel. The silica gel should be evenly stained with iodine. When the mixture is first prepared, the silica gel will appear pink, but it will turn to a rusty orange color over time.

20.8 Preparative Plates

If you use large plates (Section 20.3B), materials can be separated, and the separated components can be recovered individually from the plates. Plates used in this way are called **preparative plates**. For preparative plates, a thick layer of adsorbent is generally used. Instead of being applied as a spot or a series of spots, the mixture to be separated is applied as a line of material about 1 cm from the bottom of the plate. As the plate is developed, the separated materials form bands. After development, you can observe the separated bands, usually by UV light, and outline the zones in pencil. If the method of visualization is destructive, most of the plate is covered with paper to protect it, and the reagent is applied only at the extreme edge of the plate.

Once the zones have been identified, the adsorbent in those bands is scraped from the plate and extracted with solvent to remove the adsorbed material. Filtration removes the adsorbent, and evaporation of the solvent gives the recovered component from the mixture.

20.9 The R_f Value

Thin-layer chromatography conditions include

1. Solvent system
2. Adsorbent
3. Thickness of the adsorbent layer
4. Relative amount of material spotted

Under an established set of such conditions, a given compound always travels a fixed distance relative to the distance the solvent front travels. This ratio of the distance the compound travels to the distance the solvent travels is called the **R_f value**. The symbol R_f stands for “retardation factor,” or “ratio-to-front,” and it is expressed as a decimal fraction:


$$R_f = \frac{\text{distance traveled by substance}}{\text{distance traveled by solvent front}}$$

When the conditions of measurement are completely specified, the R_f value is constant for any given compound, and it corresponds to a physical property of that compound.

The R_f value can be used to identify an unknown compound, but like any other identification based on a single piece of data, the R_f value is best confirmed with some additional data. Many compounds can have the same R_f value, just as many compounds have the same melting point.

It is not always possible, in measuring an R_f value, to duplicate exactly the conditions of measurement another researcher has used. Therefore, R_f values tend to be of more use to a single researcher in one laboratory than they are to researchers in different laboratories. The only exception to this is when two researchers use TLC plates from the same source, as in commercial plates, or know the *exact* details of how the plates were prepared. Nevertheless, the R_f value can be a useful guide. If exact values cannot be relied on, the relative values can provide another researcher with useful information about what to expect. Anyone using published R_f values will find it a good idea to check them by comparing them with standard substances whose identity and R_f values are known.

To calculate the R_f value for a given compound, measure the distance that the compound has traveled from the point at which it was originally spotted. For spots that are not too large, measure to the center of the migrated spot. For large spots, the measurement should be repeated on a new plate, using less material. For spots that show tailing, the measurement is made to the “center of gravity” of the spot. This first distance measurement is then divided by the distance the solvent front has traveled from the same original spot. A sample calculation of the R_f values of two compounds is illustrated in Figure 20.8.

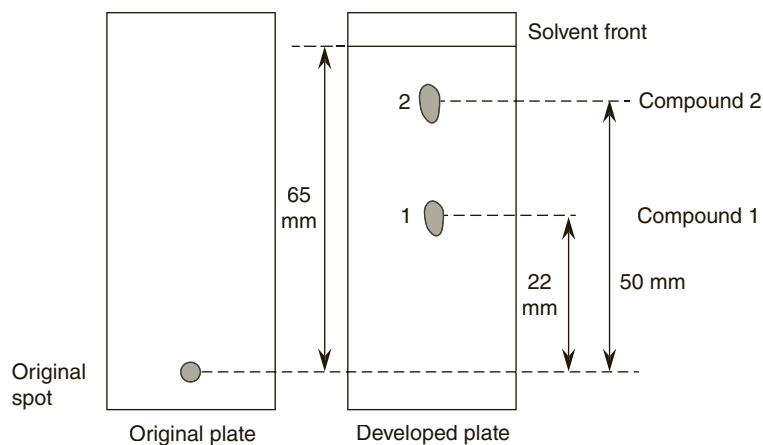

$$R_f(\text{compound 1}) = \frac{22}{65} = 0.34$$

$$R_f(\text{compound 2}) = \frac{50}{65} = 0.77$$

Figure 20.8
Sample calculation of R_f values.

20.10 Thin-Layer Chromatography Applied in Organic Chemistry

Thin-layer chromatography has several important uses in organic chemistry. It can be used in the following applications:

1. To establish that two compounds are identical
2. To determine the number of components in a mixture
3. To determine the appropriate solvent for a column-chromatographic separation
4. To monitor a column-chromatographic separation
5. To check the effectiveness of a separation achieved on a column, by crystallization or by extraction
6. To monitor the progress of a reaction

In all these applications, TLC has the advantage that only small amounts of material are necessary. Material is not wasted. With many of the visualization methods, less than a tenth of a microgram (10^{-7} g) of material can be detected. On the other hand, samples as large as a milligram may be used. With preparative plates that are large (about 9 in. on a side) and have a relatively thick coating of adsorbent ($>500 \mu\text{m}$), it is often possible to separate from 0.2 to 0.5 g of material at one time. The main disadvantage of TLC is that volatile materials cannot be used, because they evaporate from the plates.

Thin-layer chromatography can establish that two compounds suspected to be identical are in fact identical. Simply spot both compounds side by side on a single plate and develop the plate. If both compounds travel the same distance on the plate (have the same R_f value), they are probably identical. If the spot positions are not the same, the compounds are definitely not identical. It is important to spot compounds *on the same plate*. This is especially important with hand-dipped microscope slides. Because they vary widely from plate to plate, no two plates have exactly the same thickness of adsorbent. If you use commercial plates, this precaution is not necessary, although it is nevertheless a good idea.

Thin-layer chromatography can establish whether a sample is a single substance or a mixture. A single substance gives a single spot no matter what solvent is used to develop the plate. On the other hand, the number of components in a mixture can be established by trying various solvents on a mixture. A word of caution

should be given. It may be difficult, in dealing with compounds of very similar properties, isomers for example, to find a solvent that will separate the mixture. Inability to achieve a separation is not absolute proof that a sample is a single pure substance. Many compounds can be separated only by *multiple developments* of the TLC slide with a fairly nonpolar solvent. In this method, you remove the plate after the first development and allow it to dry. After being dried, it is placed in the chamber again and developed once more. This effectively doubles the length of the slide. At times, several developments may be necessary.

When a mixture is to be separated, you can use TLC to choose the best solvent to separate it if column chromatography is contemplated. You can try various solvents on a plate coated with the same adsorbent as will be used in the column. The solvent that resolves the components best will probably work well on the column. These small-scale experiments are quick, use very little material, and save time that would be wasted by attempting to separate the entire mixture on the column. Similarly, TLC plates can **monitor** a column. A hypothetical situation is shown in Figure 20.9.

A solvent was found that would separate the mixture into four components (A–D). A column was run using this solvent, and 11 fractions of 15 mL each were collected. Thin-layer analysis of the various fractions showed that Fractions 1–3 contained Component A; Fractions 4–7, Component B; Fractions 8–9, Component C; and Fractions 10–11, Component D. A small amount of cross-contamination was observed in Fractions 3, 4, 7, and 9.

In another TLC example, a researcher found a product from a reaction to be a mixture. It gave two spots, A and B, on a TLC slide. After the product was crystallized, the crystals were found by TLC to be pure A, whereas the mother liquor was found to have a mixture of A and B. The crystallization was judged to have purified A satisfactorily.

Finally, it is often possible to monitor the progress of a reaction by TLC. At various points during a reaction, samples of the reaction mixture are taken and subjected to TLC analysis. An example is given in Figure 20.10. In this case, the desired reaction was the conversion of A to B. At the beginning of the reaction (0 hr), a TLC slide was prepared that was spotted with pure A, pure B, and the reaction mixture. Similar slides were prepared at 0.5, 1, 2, and 3 hours after the start of the reaction. The slides showed that the reaction was complete in 2 hours. When the reaction was run longer than 2 hours, a new compound, side product C, began to appear. Thus, the optimum reaction time was judged to be 2 hours.

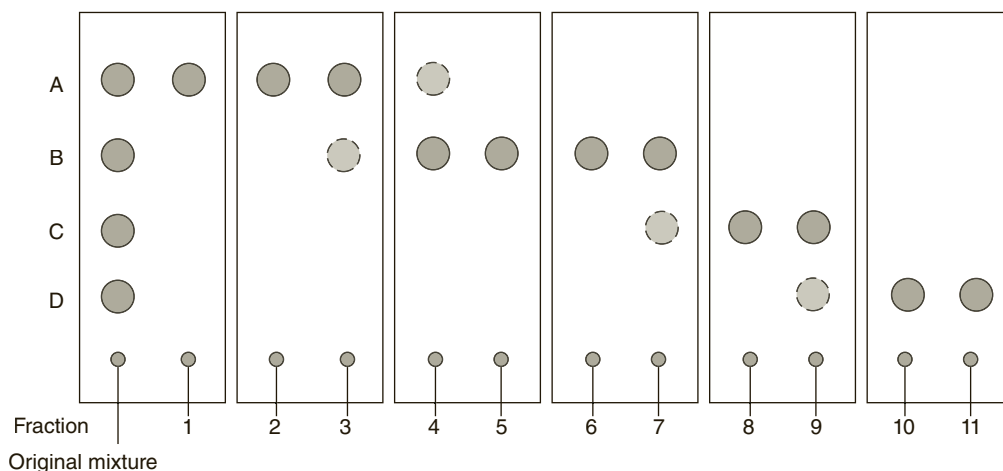


Figure 20.9
Monitoring a column.

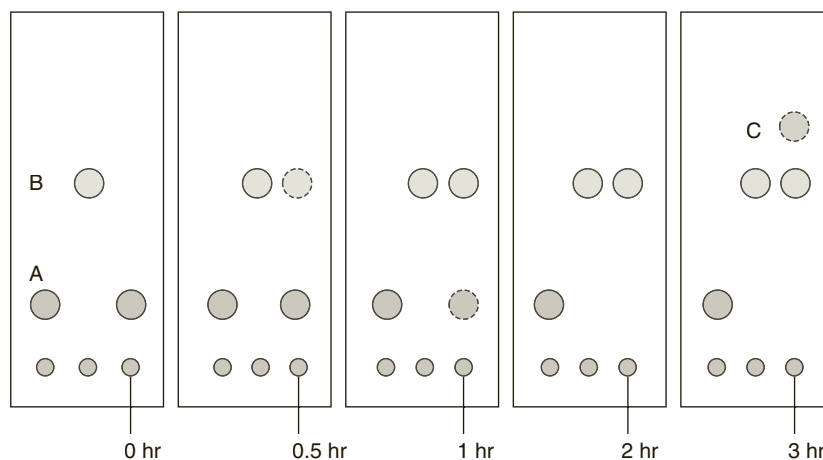


Figure 20.10
Monitoring a reaction.

20.11 Paper Chromatography

Paper chromatography is often considered to be related to thin-layer chromatography. The experimental techniques are somewhat like those of TLC, but the principles are more closely related to those of extraction. Paper chromatography is actually a liquid–liquid partitioning technique, rather than a solid–liquid technique. For paper chromatography, a spot is placed near the bottom of a piece of high-grade filter paper (Whatman No. 1 is often used). Then the paper is placed in a developing chamber. The development solvent ascends the paper by capillary action and moves the components of the spotted mixture upward at differing rates. Although paper consists mainly of pure cellulose, the cellulose itself does not function as the stationary phase. Rather, the cellulose absorbs water from the atmosphere, especially from an atmosphere saturated with water vapor. Cellulose can absorb up to about 22% of water. It is this water adsorbed on the cellulose that functions as the stationary phase. To ensure that the cellulose is kept saturated with water, many development solvents used in paper chromatography contain water as a component. As the solvent ascends the paper, the compounds are partitioned between the stationary water phase and the moving solvent. Because the water phase is stationary, the components in a mixture that are most highly water-soluble, or those that have the greatest hydrogen-bonding capacity, are the ones that are held back and move most slowly. Paper chromatography applies mostly to highly polar compounds or to those that are polyfunctional. The most common use of paper chromatography is for sugars, amino acids, and natural pigments. Because filter paper is manufactured consistently, R_f values can often be relied on in paper chromatographic work. However, R_f values are customarily measured from the leading edge (top) of the spot—not from its center, as is customary in TLC.

PROBLEMS

1. A student spots an unknown sample on a TLC plate and develops it in dichloromethane solvent. Only one spot, for which the R_f value is 0.95, is observed. Does this indicate that the unknown material is a pure compound? What can be done to verify the purity of the sample?
2. You and another student were each given an unknown compound. Both samples contained colorless material. You each used the same brand of commercially prepared TLC plate and developed the plates using the same solvent.

Each of you obtained a single spot of $R_f = 0.75$. Were the two samples necessarily the same substance? How could you prove unambiguously that they were identical using TLC?

3. Consider a sample that is a mixture composed of biphenyl, benzoic acid, and benzyl alcohol. The sample is spotted on a TLC plate and developed in a dichloromethane–cyclohexane solvent mixture. Predict the *relative* R_f values for the three components in the sample. *Hint:* See Table 19.3.
4. Calculate the R_f value of a spot that travels 5.7 cm, with a solvent front that travels 13 cm.
5. A student spots an unknown sample on a TLC plate and develops it in pentane solvent. Only one spot, for which the R_f value is 0.05, is observed. Is the unknown material a pure compound? What can be done to verify the purity of the sample?
6. A *colorless* unknown substance is spotted on a TLC plate and developed in the correct solvent. The spots do not appear when visualization with a UV lamp or iodine vapors is attempted. What could you do in order to visualize the spots if the compound is
 - a. An alkyl halide
 - b. A ketone
 - c. An amino acid
 - d. A sugar

High-Performance Liquid Chromatography (HPLC)

The separation that can be achieved is greater if the column packing used in column chromatography is made more dense by using an adsorbent that has a smaller particle size. The solute molecules encounter a much larger surface area on which they can be adsorbed as they pass through the column packing. At the same time, the solvent spaces between the particles are reduced in size. As a result of this tight packing, equilibrium between the liquid and solid phases can be established rapidly with a fairly short column, and the degree of separation is markedly improved. The disadvantage of making the column packing more dense is that the solvent flow rate becomes very slow or even stops. Gravity is not strong enough to pull the solvent through a tightly packed column.

A recently developed technique can be applied to obtain much better separations with tightly packed columns. A pump forces the solvent through the column packing. As a result, solvent flow rate is increased and the advantage of better separation is retained. This technique, called **high-performance liquid chromatography (HPLC)**, is becoming widely applied to problems where separations by ordinary column chromatography are unsatisfactory. Because the pump often provides pressures in excess of 1000 pounds per square inch (psi), this method is also known as **high-pressure liquid chromatography**. High pressures are not required, however, and satisfactory separations can be achieved with pressures as low as 100 psi.

The basic design of an HPLC instrument is shown in Figure 21.1. The instrument contains the following essential components:

1. Solvent reservoir
2. Solvent filter and degasser
3. Pump
4. Pressure gauge
5. Sample injection system
6. Column
7. Detector
8. Amplifier and electronic controls
9. Chart recorder

There may be other variations on this simple design. Some instruments have heated ovens in order to maintain the column at a specified temperature, fraction collectors, and microprocessor-controlled data-handling systems. Additional filters for the solvent and sample may also be included. You may find it interesting to compare this schematic diagram with that shown in Technique 22 (Figure 22.2) for a gas-chromatography instrument. Many of the essential components are common to both types of instruments.

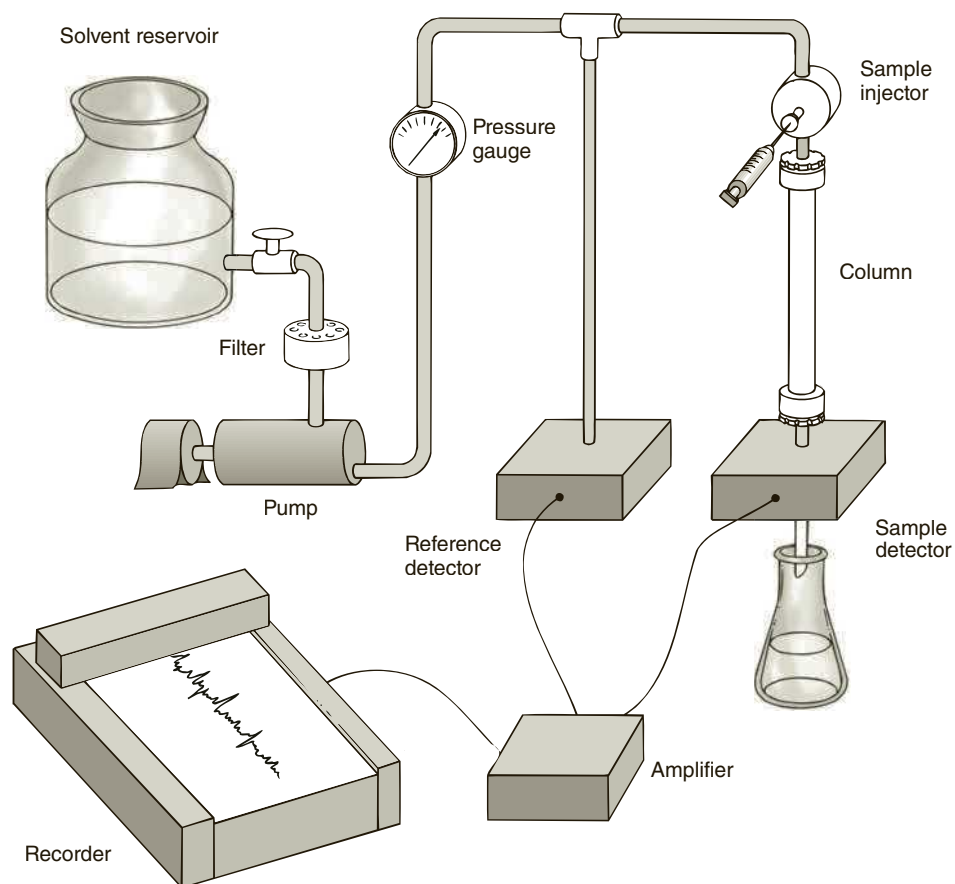


Figure 21.1
Schematic diagram of a high-performance liquid chromatograph.

21.1 Adsorbents and Columns

The most important factor to consider when choosing a set of experimental conditions is the nature of the material packed into the column. You must also consider the size of the column that will be selected. The chromatography column is generally packed with silica or alumina adsorbents. Unlike column chromatography, however, the adsorbents used for HPLC have a much smaller particle size. Typically, particle size ranges from 5 to 20 μm in diameter for HPLC; it is on the order of 100 μm for column chromatography.

The adsorbent is packed into a column that can withstand the elevated pressures typical of this type of experiment. Generally, the column is constructed of stainless steel, although some columns that are constructed of a rigid polymeric material ("PEEK" —Poly Ether Ether Ketone) are available commercially. A strong column is required to withstand the high pressures that may be used. The columns are fitted with stainless steel connectors, which ensure a pressure-tight fit between the tubing that connects the column to the other components of the instrument.

Columns that fulfill a large number of specialized purposes are available. In this chapter, we consider only the four important types of columns. These are

1. Normal-phase chromatography
2. Reversed-phase chromatography
3. Ion-exchange chromatography
4. Size-exclusion chromatography

In most types of chromatography, the adsorbent is more polar than the mobile phase. For example, the solid packing material, which may be either silica or alumina, has a stronger affinity for polar molecules than does the solvent. As a result, the molecules in the sample adhere strongly to the solid phase, and their progress down the column is much slower than the rate at which solvent moves through the column. The time required for a substance to move through the column can be altered by changing the polarity of the solvent. In general, as the solvent becomes more polar, the faster substances move through the column. This type of behavior is known as **normal-phase chromatography**. In HPLC, you inject a sample onto a normal-phase column and elute it by varying the polarity of the solvent, much as you do with ordinary column chromatography. Disadvantages of normal-phase chromatography are that retention times tend to be long, and bands have a tendency to “tail.”

These disadvantages can be ameliorated by selecting a column in which the solid support is *less polar* than the moving solvent phase. This type of chromatography is known as **reversed-phase chromatography**. In this type of chromatography, the silica column packing is treated with alkylating agents. As a result, nonpolar alkyl groups are bonded to the silica surface, making the adsorbent nonpolar. The alkylating agents that are used most commonly can attach methyl ($-\text{CH}_3$), octyl ($-\text{C}_8\text{H}_{17}$), or octadecyl ($-\text{C}_{18}\text{H}_{37}$) groups to the silica surface. The latter variation, where an 18-carbon chain is attached to the silica, is the most popular. This type of column is known as a **C-18 column**. The bonded alkyl groups have an effect similar to that which would be produced by an extremely thin organic solvent layer coating the surface of the silica particles. The interactions that take place between the substances dissolved in the solvent and the stationary phase thus become more like those observed in a liquid-liquid extraction. The solute particles distribute themselves between the two “solvents”—that is, between the moving solvent and the organic coating on the silica. The longer the chains of the alkyl groups that are bonded to the silica, the more effective the alkyl groups are as they interact with solute molecules.

Reversed-phase chromatography is widely used because the rate at which solute molecules exchange between moving phase and stationary phase is rapid, which means that substances pass through the column relatively quickly. Furthermore, problems arising from the “tailing” of peaks are reduced. A disadvantage of this type of column, however, is that the chemically bonded solid phases tend to decompose. The organic groups are slowly hydrolyzed from the surface of the silica, which leaves a normal silica surface exposed. Thus, the chromatographic process that takes place on the column slowly shifts from a reversed-phase to a normal-phase separation mechanism.

Another type of solid support that is sometimes used in reversed-phase chromatography is organic polymer beads. These beads present a surface to the moving phase that is largely organic in nature.

For solutions of ions, select a column that is packed with an ion-exchange resin. This type of chromatography is known as **ion-exchange chromatography**. The ion-exchange resin that is chosen can be either an anion-exchange resin or a cation-exchange resin, depending upon the nature of the sample being examined.

A fourth type of column is known as a **size-exclusion column** or a **gel-filtration column**. The interaction that takes place on this type of column is similar to that described in Technique 19, Section 19.15.

21.2 Column Dimensions

The dimensions of the column that you use depend upon the application. For analytical applications, a typical column is constructed of tubing that has an inside diameter of between 4 and 5 mm, although analytical columns with inside diameters

of 1 or 2 mm are also available. A typical analytical column has a length of about 7.5 to 30 cm. This type of column is suitable for the separation of a 0.1- to 5-mg sample. With columns of smaller diameter, it is possible to perform an analysis with samples smaller than 1 *microgram*.

High-performance liquid chromatography is an excellent analytical technique, but the separated compounds may also be isolated. The technique can be used for preparative experiments. Just as in column chromatography, the fractions can be collected into individual receiving containers as they pass through the column. The solvents can be evaporated from these fractions, allowing you to isolate separated components of the original mixture. Samples that range in size from 5 to 100 mg can be separated on a semipreparative, or **semiprep**, column. The dimensions of a semiprep column are typically 8 mm inside diameter and 10 cm in length. A semiprep column is a practical choice when you wish to use the same column for both analytical and preparative separations. A semiprep column is small enough to provide reasonable sensitivity in analyses, but it is also capable of handling moderate-size samples when you need to isolate the components of a mixture. Even larger samples can be separated using a **preparative column**. This type of column is useful when you wish to collect the components of a mixture and then use the pure samples for additional study (e.g., for a subsequent chemical reaction or for spectroscopic analysis). A preparative column may be as large as 20 mm in inside diameter and 30 cm in length. A preparative column can handle samples as large as 1 g per injection.

21.3 Solvents

The choice of solvent used for an HPLC separation depends on the type of chromatographic process selected. For a normal-phase separation, the solvent is selected based on its polarity. The criteria described in Technique 19, Section 19.4B, are used. A solvent of very low polarity might be pentane, petroleum ether, hexane, or carbon tetrachloride; a solvent of very high polarity might be water, acetic acid, methanol, or 1-propanol. For a reversed-phase experiment, a less polar solvent causes solutes to migrate *faster*. For example, for a mixed methanol-water solvent, as the percentage of methanol in the solvent increases (solvent becomes less polar), the time required to elute the components of a mixture from a column decreases. The behavior of solvents as eluents in a reversed-phase chromatography would be the reverse of the order shown in Table 19.2 (Technique 19).

If a single solvent (or solvent mixture) is used for the entire separation, the chromatogram is said to be **isochratic**. Special electronic devices are available with HPLC instruments that allow you to program changes in the solvent composition from the beginning to the end of the chromatography. These are called **gradient elution systems**. With gradient elution, the time required for a separation may be shortened considerably.

The need for pure solvents is especially acute with HPLC. The narrow bore of the column and the very small particle size of the column packing require that solvents be particularly pure and free of insoluble residue. In most cases, the solvents must be filtered through ultrafine filters and **degassed** (have dissolved gases removed) before they can be used.

The solvent gradient is chosen so that the eluting power of the solvent increases over the duration of the experiment. The result is that components of the mixture that tend to move very slowly through the column are caused to move faster as the eluting power of the solvent gradually increases. The instrument can be programmed to change the composition of the solvent following a linear gradient or a nonlinear gradient, depending on the specific requirements of the separation.

21.4 Detectors

A flow-through **detector** must be provided to determine when a substance has passed through the column. In most applications, the detector detects either the change in index of refraction of the liquid as its composition changes or the presence of solute by its absorption of ultraviolet or visible light. The signal generated by the detector is amplified and treated electronically in a manner similar to that found in gas chromatography (Technique 22, Section 22.6).

A detector that responds to changes in the index of refraction of the solution may be considered the most universal of the HPLC detectors. The refractive index of the liquid passing through the detector changes slightly, but significantly, as the liquid changes from pure solvent to a liquid where the solvent contains some type of organic solute. This change in refractive index can be detected and compared to the refractive index of pure solvent. The difference in index values is then recorded as a peak on a chart. A disadvantage of this type of detector is that it must respond to very small changes in refractive index. As a result, the detector tends to be unstable and difficult to balance.

When the components of the mixture have some type of absorption in the ultraviolet or visible regions of the spectrum, a detector that is adjusted to detect absorption at a particular wavelength of light can be used. This type of detector is much more stable, and the readings tend to be more reliable. Unfortunately, many organic compounds do not absorb ultraviolet light, and this type of detector cannot be used.

21.5 Presentation of Data

The data produced by an HPLC instrument appear in the form of a chart, where detector response is the vertical axis and time is represented on the horizontal axis. These are recorded on a continuously moving strip of chart paper, although they may also be observed in graphic form on a computer display. In virtually all respects, the form of the data is identical to that produced by a gas chromatograph; in fact, in many cases, the data-handling system for the two types of instruments is essentially identical. To understand how to analyze the data from an HPLC instrument, read Sections 22.12 and 22.13 in Technique 22.

REFERENCE

Rubinson, K. A. *Chemical Analysis*. Little, Brown and Co: Boston, 1987. Chapter 14, Liquid Chromatography.

PROBLEMS

1. For a mixture of biphenyl, benzoic acid, and benzyl alcohol, predict the order of elution and describe any differences that you would expect for a normal-phase HPLC experiment (in hexane solvent) compared with a reversed-phase experiment (in tetrahydrofuran-water solvent).
2. How would the *gradient elution program* differ between normal-phase and reversed-phase chromatography?

Gas Chromatography

Gas chromatography is one of the most useful instrumental tools for separating and analyzing organic compounds that can be vaporized without decomposition. Common uses include testing the purity of a substance and separating the components of a mixture. The relative amounts of the components in a mixture may also be determined. In some cases, gas chromatography can be used to identify a compound. In microscale work, it can also be used as a preparative method to isolate pure compounds from a small amount of a mixture.

Gas chromatography resembles column chromatography in principle, but it differs in three respects. First, the partitioning processes for the compounds to be separated are carried out between a **moving gas phase** and a **stationary liquid phase**. (Recall that in column chromatography the moving phase is a liquid and the stationary phase is a solid adsorbent.) A second difference is that the temperature of the gas system can be controlled because the column is contained in an insulated oven. And third, the concentration of any given compound in the gas phase is a function of its vapor pressure only. Because gas chromatography separates the components of a mixture primarily on the basis of their vapor pressures (or boiling points), this technique is also similar in principle to fractional distillation. In microscale work, it is sometimes used to separate and isolate compounds from a mixture; fractional distillation would normally be used with larger amounts of material.

Gas chromatography (GC) is also known as **vapor-phase chromatography** (VPC) and as **gas-liquid partition chromatography** (GLPC). All three names, as well as their indicated abbreviations, are often found in the literature of organic chemistry. In reference to the technique, the last term, GLPC, is the most strictly correct and is preferred by most authors.

22.1 The Gas Chromatograph

The apparatus used to carry out a gas-liquid chromatographic separation is generally called a **gas chromatograph**. A typical student-model gas chromatograph, the GOW-MAC model 69-350, is illustrated in Figure 22.1. A schematic block diagram of a basic gas chromatograph is shown in Figure 22.2. The basic elements of the apparatus are apparent. In short, the sample is injected into the chromatograph, and it is immediately vaporized in a heated injection chamber and introduced into a moving stream of gas, called the **carrier gas**. The vaporized sample is then swept into a column filled with particles coated with a liquid adsorbent. The column is contained in a temperature-controlled oven. As the sample passes through the column, it is subjected to many gas-liquid partitioning processes, and the components are separated. As each component leaves the column, its presence is detected by an electrical detector that generates a signal that is recorded on a strip chart recorder.

Many modern instruments are also equipped with a microprocessor, which can be programmed to change parameters, such as the temperature of the oven, while a mixture is being separated on a column. With this capability, it is possible to optimize the separation of components and to complete a run in a relatively short time.

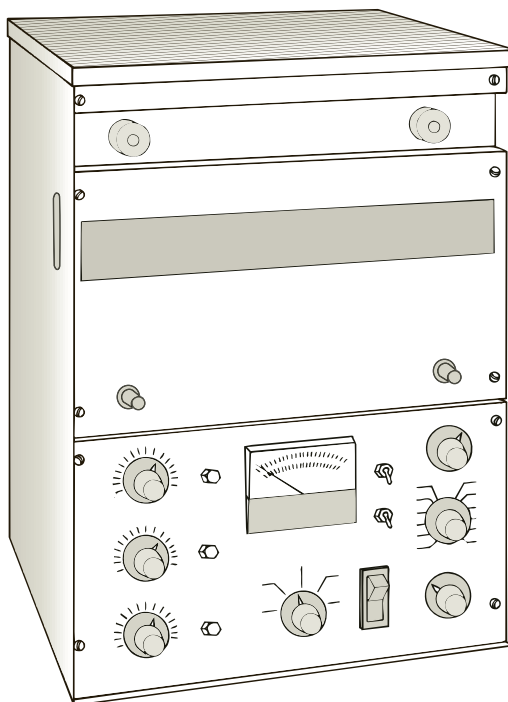


Figure 22.1
Gas chromatograph.

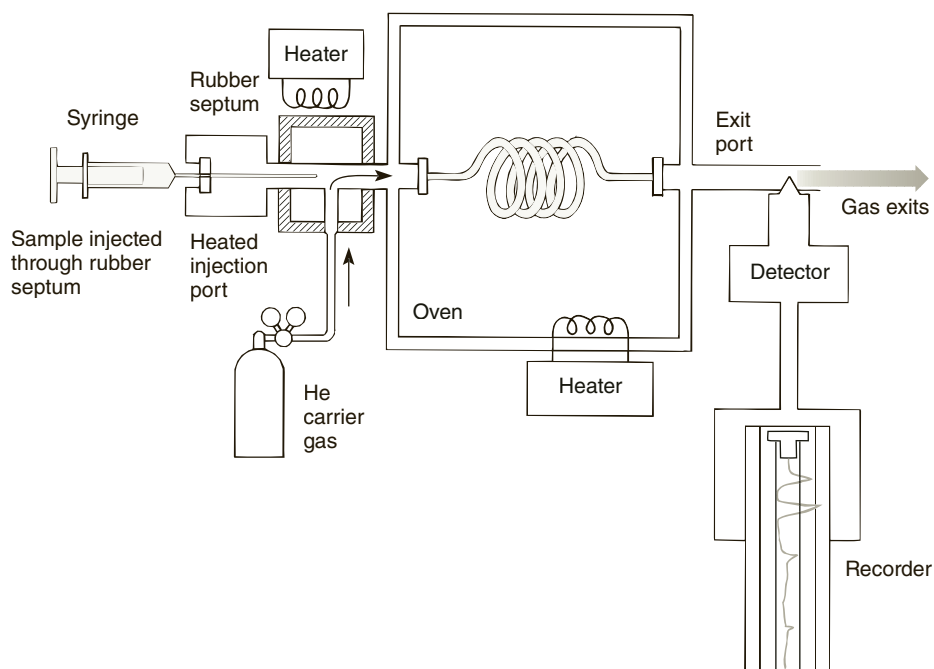


Figure 22.2
Schematic diagram of a gas chromatograph.

22.2 The Column

The heart of the gas chromatograph is the column. There are two types in common use: packed and capillary columns.

Packed Columns. These columns are usually constructed of stainless steel tubing with diameters of 1/8 in. (3 mm) or 1/4 in. (6 mm) and lengths of from 4 to 12 feet. The column is packed with a liquid or low-melting solid as the **stationary phase** distributed on a solid **support material**. The stationary phase must be relatively nonvolatile, that is, it should have a low vapor pressure and a high boiling point. Some typical stationary phases used with packed columns are listed in Table 22.1. Typical support materials are shown in Table 22.2. The most common support material consists of diatomaceous earth (Chromosorb).

Packed columns are bought from commercial sources or may sometimes be made in the laboratory by researchers. Basically, you dissolve one of the stationary phases listed in Table 22.1 in methylene chloride. Then you add the support material to the solution followed by removal of the solvent on a rotary evaporator (see Technique 7, Section 7.11, and Figure 7.19). The evaporation process evenly distributes the stationary phase onto the support material and yields a dry solid. In the final step, the solid, consisting of the stationary phase coated on the support

Table 22.1 Typical Stationary Phases

	Type	Composition	Maximum Temperature (°C)	Typical Use
	Apiezon (L, M, N, etc.)	Hydrocarbon greases (varying MW)	250–300	Hydrocarbons
Increasing polarity ↓	SE-30	Methyl silicone rubber	350	General applications
	DC-200	Silicone oil (R=CH ₃)	225	Aldehydes, ketones, halocarbons
		$\text{R}_3\text{Si}-\text{O}-\left[\begin{array}{c} \text{R} \\ \\ \text{Si}-\text{O} \\ \\ \text{R} \end{array} \right]_n-\text{SiR}_3$		
	DC-710	Silicone oil (R=CH ₃) (R=C ₆ H ₅)	300	General applications
		$\left[\begin{array}{c} \text{R}' \\ \\ \text{Si}-\text{O} \\ \\ \text{R} \end{array} \right]_n-$		
Carbowaxes (400–20M)	Polyethylene glycols (varying chain lengths)	Polyether HO—(CH ₂ CH ₂ —O) _n —CH ₂ CH ₂ OH	Up to 250	Alcohols, ethers, halocarbons
DEGS	Diethylene glycol succinate	Polyester $\left(\text{CH}_2\text{CH}_2-\text{O}-\overset{\text{O}}{\parallel}{\text{C}}-(\text{CH}_2)_2-\overset{\text{O}}{\parallel}{\text{C}}-\text{O} \right)_n$	200	General applications

Table 22.2 Typical Solid Supports

Crushed firebrick	Chromosorb T
Nylon beads	(Teflon beads)
Glass beads	Chromosorb P
Silica	(pink diatomaceous earth,
Alumina	high absorptivity, pH 6–7)
Charcoal	Chromosorb W
Molecular sieves	(white diatomaceous earth,
	medium absorptivity, pH 8–10)
	Chromosorb G
	(like the above,
	low absorptivity, pH 8.5)

material, is packed into the stainless steel tubing as evenly as possible. After plugging the ends of the tubing with glass wool to prevent the solid from coming out, the tubing is rolled into a coil that will fit into the oven of the gas chromatograph with the two ends connected to the gas entrance and exit ports (see Figure 22.2).

The selection of the packed column depends on the application. If one wants to separate nonpolar compounds that vary only by boiling point, one often uses one of the polydimethyl siloxanes (methyl silicones) columns such as SE-30 or DC-200. For more polar compounds, chemists will select a silicone column that has attached methyl and phenyl groups on the silicone polymer (DC-710). Separating even more polar compounds calls for a polyethylene glycol (Carbowax) column or a column packed with diethylene glycol succinate (DEGS). Chemists will need to carefully note the maximum temperature that can be employed with the columns (Table 22.1). Above the specified temperature, the liquid phase itself will begin to “bleed” off the column.

Packed columns are cheaper to buy and can separate larger quantities of material than capillary columns. Packed columns, however, are not as efficient in separating materials, especially with compounds with very similar polarities or boiling points.

Capillary Columns. Many gas chromatographs sold today use capillary columns rather than packed columns. Capillary columns are made of very thin fused-silica with an inner diameter of about 0.25 mm. Typically the columns are very long, often 25 meters up to 100 meters in length. The stationary phase is coated as a film on the inner surface with a thickness of about 0.25 μm . The film is bonded to the silica and cross-linked to improve thermal stability and to help prevent bleeding of the stationary phase from the column. Most of the capillary columns do not have any support material in them.

The liquid stationary phases used with capillary columns are similar to those used in packed columns. The most common stationary phases are polysiloxanes (silicones) that contain various substituents that modify the polarity of the phase. Polydimethyl siloxane (methyl silicone) is nonpolar. Replacing methyl groups with increasing numbers of phenyl substituents increases the polarity of the silicone. For example, J&W DB-1, or a similar stationary phase sold by other companies,¹ has the same properties as methyl silicone (SE-30 or DC-200) and is used with nonpolar

¹ Many companies supply capillary columns: J&W (Alltech), Supelco, HP, Chrompack, and Quadrex are some of the suppliers. Comparison of chemical compositions, polarities, and applications may be made by consulting: <http://www.quadrexcorp.com/2009/index.htm> (accessed April 23, 2011).

compounds that separate by boiling point differences. For other applications, a J&W DB-5 may be used. This silicone consists of a material with 95% methyl and 5% phenyl substituents on the silicone. A DB-17 column has 50% phenyl groups replacing methyl groups, and would have a medium polarity (similar to DC-710). J&W DB-wax is used to separate much more polar compounds and is similar to diethylene glycols (Carbowax). Chiral capillary columns are also available that can separate enantiomers (Section 22.8).

Because of the length and small diameter of capillary columns, there is increased interaction between the compounds in the mixture and the stationary phase. Capillary columns are, therefore, much more efficient in separating compounds with similar properties than with packed columns, but only if a dilute solution of a mixture of compounds is injected into the column. In order to obtain a satisfactory separation, you must dissolve about 1 drop of a mixture in about 2 mL of a solvent such as methylene chloride or pentane. *About 1 μL of this dilute sample is injected onto the column. In contrast, 1 to 10 μL of an undiluted sample can often be analyzed on a packed column without overloading the column. These expensive capillary columns must be purchased from commercial companies specializing in their manufacture. Sensible researchers would never try to make their own capillary columns!*

In the final step, the liquid-phase-coated support material is packed into the tubing as evenly as possible. The tubing is bent or coiled so that it fits into the oven of the gas chromatograph with its two ends connected to the gas entrance and exit ports.

Selection of a liquid phase usually revolves about two factors. First, most of them have an upper temperature limit above which they cannot be used. Above the specified limit of temperature, the liquid phase itself will begin to “bleed” off the column. Second, the materials to be separated must be considered. For polar samples, it is usually best to use a polar liquid phase; for nonpolar samples, a nonpolar liquid phase is indicated. The liquid phase performs best when the substances to be separated *dissolve* in it.

Most researchers today buy packed columns from commercial sources, rather than pack their own. A wide variety of types and lengths is available. Alternatives to packed columns are Golay or glass capillary columns of diameters 0.1–0.2 mm. With these columns, no solid support is required, and the liquid is coated directly on the inner walls of the tubing. Liquid phases commonly used in glass capillary columns are similar in composition to those used in packed columns. They include DB-1 (similar to SE-30), DB-17 (similar to DC-710), and DB-WAX (similar to Carbowax 20M). The length of a capillary column is usually very long, typically 50–100 ft. Because of the length and small diameter, there is increased interaction between the sample and the stationary phase. Gas chromatographs equipped with these small-diameter columns are able to separate components more effectively than instruments using larger packed columns.

22.3 Principles of Separation

After a column is selected and installed, the carrier gas (usually helium, argon, or nitrogen) is allowed to flow through the column supporting the liquid phase. The mixture of compounds to be separated is introduced into the carrier gas stream, where its components are equilibrated (or partitioned) between the moving gas phase and the stationary liquid phase (Figure 22.3). The latter is held stationary because it is adsorbed onto the surfaces of the support material.

The sample is introduced into the gas chromatograph by a microliter syringe. It is injected as a liquid or as a solution through a rubber septum into a heated chamber, called the **injection port**, where it is vaporized and mixed with the carrier gas.

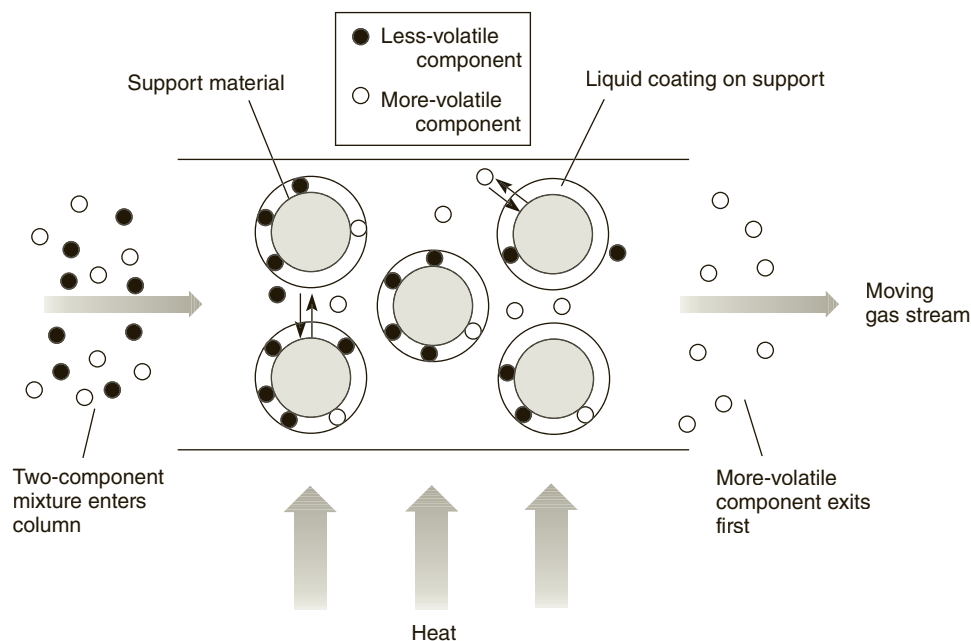


Figure 22.3
The separation process.

As this mixture reaches the column, which is heated in a controlled oven, it begins to equilibrate between the liquid and gas phases. The length of time required for a sample to move through the column is a function of how much time it spends in the vapor phase and how much time it spends in the liquid phase. The more time it spends in the vapor phase, the faster it gets to the end of the column. In most separations, the components of a sample have similar solubilities in the liquid phase. Therefore, the time the different compounds spend in the vapor phase is primarily a function of the vapor pressure, and the more volatile component arrives at the end of the column first, as illustrated in Figure 22.3. By selecting the correct temperature of the oven and the correct liquid phase, the compounds in the injected mixture travel through the column at different rates and are separated.

22.4 Factors Affecting Separation

Several factors determine the rate at which a given compound travels through a gas chromatograph. First of all, compounds with low boiling points will generally travel through the gas chromatograph faster than compounds of higher boiling points. This is because the column is heated, and low-boiling compounds always have higher vapor pressures than compounds of higher boiling point. In general, therefore, for compounds with the same functional group, the higher the molecular weight, the longer the retention time. For most molecules, the boiling point increases as the molecular weight increases. If the column is heated to a temperature that is too high, however, the entire mixture to be separated is flushed through the column at the same rate as the carrier gas, and no equilibration takes place with the liquid phase. On the other hand, at too low a temperature, the mixture dissolves in the liquid phase and never reevaporizes. Thus, it is retained on the column.

The second factor, the rate of flow of the carrier gas, is important. The carrier gas must not move so rapidly that molecules of the sample in the vapor phase cannot equilibrate with those dissolved in the liquid phase. This may result in poor separation between components in the injected mixture. If the rate of flow is too slow, however, the bands broaden significantly, leading to poor resolution (see Section 22.8).

The third factor is the choice of liquid phase used in the column. The molecular weights, functional groups, and polarities of the component molecules in the mixture to be separated must be considered when a liquid phase is being chosen. One generally uses a different type of material for hydrocarbons, for instance, than for esters. The materials to be separated should *dissolve* in the liquid. The useful temperature limit of the liquid phase selected must also be considered.

The fourth factor, the length of the column, is also important. Compounds that resemble one another closely, in general, require longer columns than dissimilar compounds. Many kinds of isomeric mixtures fit into the “difficult” category. The components of isomeric mixtures are so much alike that they travel through the column at very similar rates. You need a longer column, therefore, to take advantage of any differences that may exist.

22.5 Advantages of Gas Chromatography

All factors that have been mentioned must be adjusted by the chemist for any mixture to be separated. Considerable preliminary investigation is often required before a mixture can be separated successfully into its components by gas chromatography. Nevertheless, the advantages of the technique are many.

First, many mixtures can be separated by this technique when no other method is adequate. Second, as little as 1–10 μL ($1 \mu\text{L} = 10^{-6}\text{L}$) of a mixture can be separated by this technique. This advantage is particularly important when working at the microscale level. Third, when gas chromatography is coupled with an electronic recording device (see following discussion), the amount of each component present in the separated mixture can be estimated quantitatively.

The range of compounds that can be separated by gas chromatography extends from gases, such as oxygen (bp – 183°C) and nitrogen (bp – 196°C), to organic compounds with boiling points over 400°C. The only requirement for the compounds to be separated is that they have an appreciable vapor pressure at a temperature at which they can be separated and that they be thermally stable at this temperature.

22.6 Monitoring the Column (The Detector)

To follow the separation of the mixture injected into the gas chromatograph, it is necessary to use an electrical device called a **detector**. Two types of detectors in common use are the **thermal conductivity detector (TCD)** and the **flame-ionization detector (FID)**.

The thermal conductivity detector is simply a hot wire placed in the gas stream at the column exit. The wire is heated by constant electrical voltage. When a steady stream of carrier gas passes over this wire, the rate at which it loses heat and its electrical conductance have constant values. When the composition of the vapor stream changes, the rate of heat flow from the wire, and hence its resistance, changes. Helium, which has a higher thermal conductivity than most organic substances, is a common carrier gas. Thus, when a substance elutes in the vapor stream, the thermal conductivity of the moving gases will be lower than with helium alone. The wire then heats up, and its resistance decreases.

A typical TCD operates by difference. Two detectors are used: one exposed to the actual effluent gas and the other exposed to a reference flow of carrier gas only. To achieve this situation, a portion of the carrier gas stream is diverted *before* it enters the injection port. The diverted gas is routed through a reference column into which no sample has been admitted. The detectors mounted in the sample and reference columns are arranged so as to form the arms of a Wheatstone bridge circuit, as shown in Figure 22.4. As long as the carrier gas alone flows over both detectors, the circuit is in balance. However, when a sample elutes from the sample column, the bridge circuit becomes unbalanced, creating an electrical signal. This signal can

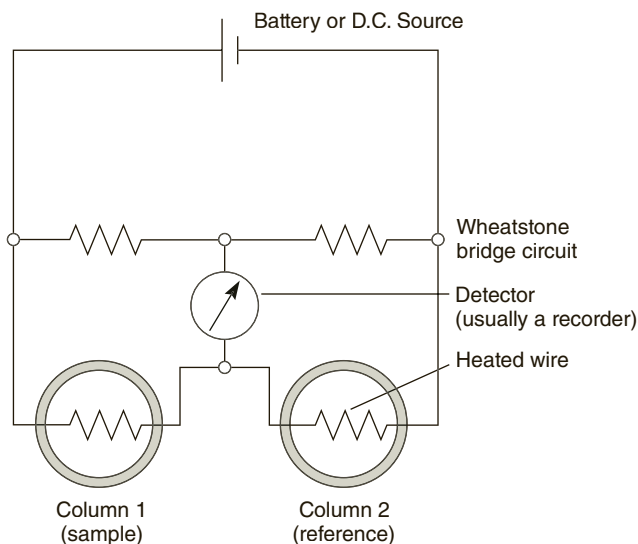


Figure 22.4
Typical thermal conductivity detector.

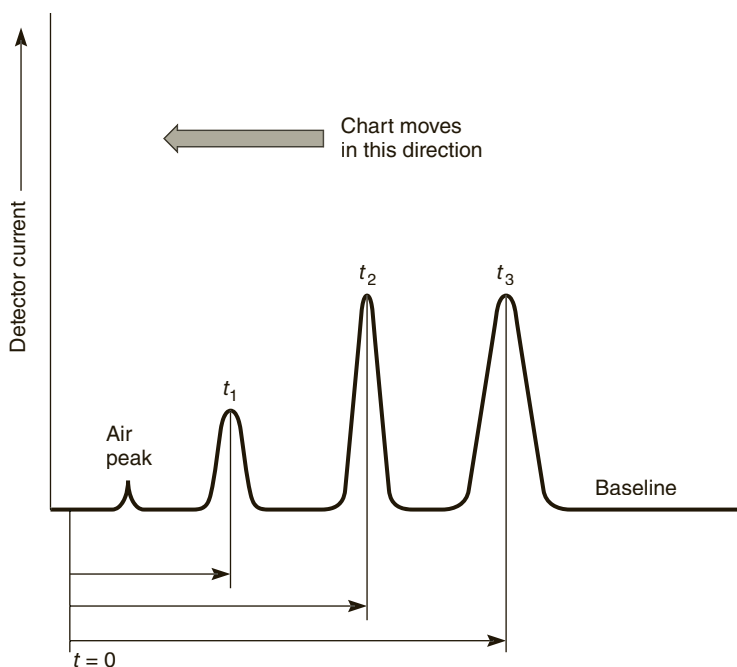


Figure 22.5
Typical chromatogram.

be amplified and used to activate a strip chart recorder. The recorder is an instrument that plots, by means of a moving pen, the unbalanced bridge current versus time on a continuously moving roll of chart paper. This record of detector response (current) versus time is called a **chromatogram**. A typical gas chromatogram is illustrated in Figure 22.5. Deflections of the pen are called **peaks**.

When a sample is injected, some air (CO_2 , H_2O , N_2 , and O_2) is introduced along with the sample. The air travels through the column almost as rapidly as the carrier gas; as it passes the detector, it causes a small pen response, thereby giving a peak,

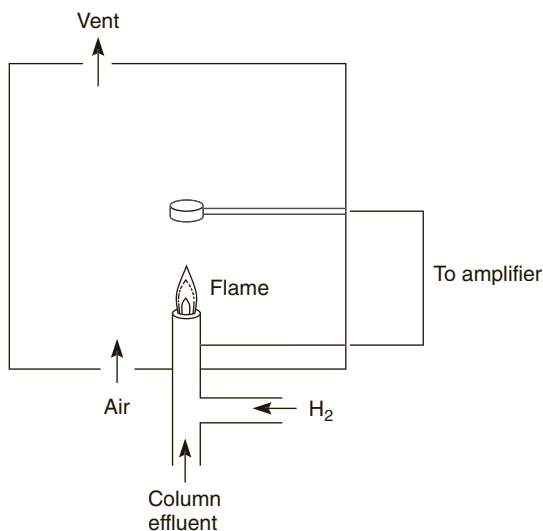


Figure 22.6
Flame-ionization detector.

called the **air peak**. At later times (t_1, t_2, t_3), the components also give rise to peaks on the chromatogram as they pass out of the column and past the detector.

In a flame-ionization detector, the effluent from the column is directed into a flame produced by the combustion of hydrogen, as illustrated in Figure 22.6. As organic compounds burn in the flame, ion fragments are produced that collect on the ring above the flame. The resulting electrical signal is amplified and sent to a recorder in a similar manner to a TCD, except that an FID does not produce an air peak. The main advantage of the FID is that it is more sensitive and can be used to analyze smaller quantities of sample. Also, because an FID does not respond to water, a gas chromatograph with this detector can be used to analyze aqueous solutions. Two disadvantages are that it is more difficult to operate and the detection process destroys the sample. Therefore, an FID gas chromatograph cannot be used to do preparative work, which is often desired in the microscale laboratory.

22.7 Retention Time

The period following injection that is required for a compound to pass through the column is called the **retention time** of that compound. For a given set of constant conditions (flow rate of carrier gas, column temperature, column length, liquid phase, injection port temperature, carrier), the retention time of any compound is always constant (much like the R_f value in thin-layer chromatography, as described in Technique 20, Section 20.9). The retention time is measured from the time of injection to the time of maximum pen deflection (detector current) for the component being observed. This value, when obtained under controlled conditions, can identify a compound by a direct comparison of it with values for known compounds determined under the same conditions. For easier measurement of retention times, most strip chart recorders are adjusted to move the paper at a rate that corresponds to time divisions calibrated on the chart paper. The retention times (t_1, t_2, t_3) are indicated in Figure 22.5 for the three peaks illustrated.

Most modern gas chromatographs are attached to a “data station,” which uses a computer or a microprocessor to process the data. With these instruments, the chart often does not have divisions. Instead, the computer prints the retention time, usually to the nearest 0.01 minute, above each peak.

22.8 Chiral Stationary Phases

A recent innovation in gas chromatography is to use chiral adsorbent materials to achieve separations of stereoisomers. The interaction between a particular stereoisomer and the chiral adsorbent may be different from the interaction between the opposite stereoisomer and the same chiral adsorbent. As a result, retention times for the two stereoisomers are likely to be sufficiently different to allow for a clean separation. The interactions between a chiral substance and the chiral adsorbent will include hydrogen-bonding and dipole-dipole attraction forces, although other properties may also be involved. One enantiomer should interact more strongly with the adsorbent than its opposite form. Thus, one enantiomer should pass through the gas chromatography column more slowly than its opposite form.

The ability of chiral adsorbents to separate stereoisomers is rapidly finding many useful applications, particularly in the synthesis of pharmaceutical agents. The biological activity of chiral substances often depends upon their stereochemistry because the living body is a highly chiral environment. A large number of pharmaceutical compounds have two enantiomeric forms that in many cases show significant

differences in their behavior and activity. The ability to prepare enantiomerically pure drugs is very important because these pure substances are much more potent (and often have fewer side effects) than their racemic analogues.

Another type of stationary phase in gas chromatography is based on molecules such as the **cyclodextrins**. With these materials, the discrimination between enantiomers depends on the interactions between the stereoisomers and the chiral cavity that is formed within these materials. Because enantiomers differ in shape, they will fit differently within the chiral cavity. The result will be that the enantiomers will pass through the cyclodextrin stationary phase at different rates, thus leading to a separation.

The cyclodextrins owe their specificity to their structure, which is based on polymers of D-(+)-glucose. The hydroxyl groups of the glucose have been alkylated, so that the cavity is relatively nonpolar. The exterior hydroxyl groups of the cyclodextrins have also been substituted with *tert*-butyldimethylsilyl groups. The result is a material that can also utilize differences in hydrogen-bonding and dipole—dipole interactions to separate stereoisomers.

The structure of one important cyclodextrin-based chiral adsorbent is shown in Figure 22.7. Gas chromatography using this chiral adsorbent as a stationary phase has been used to separate a wide variety of stereoisomers. In one recent publication,


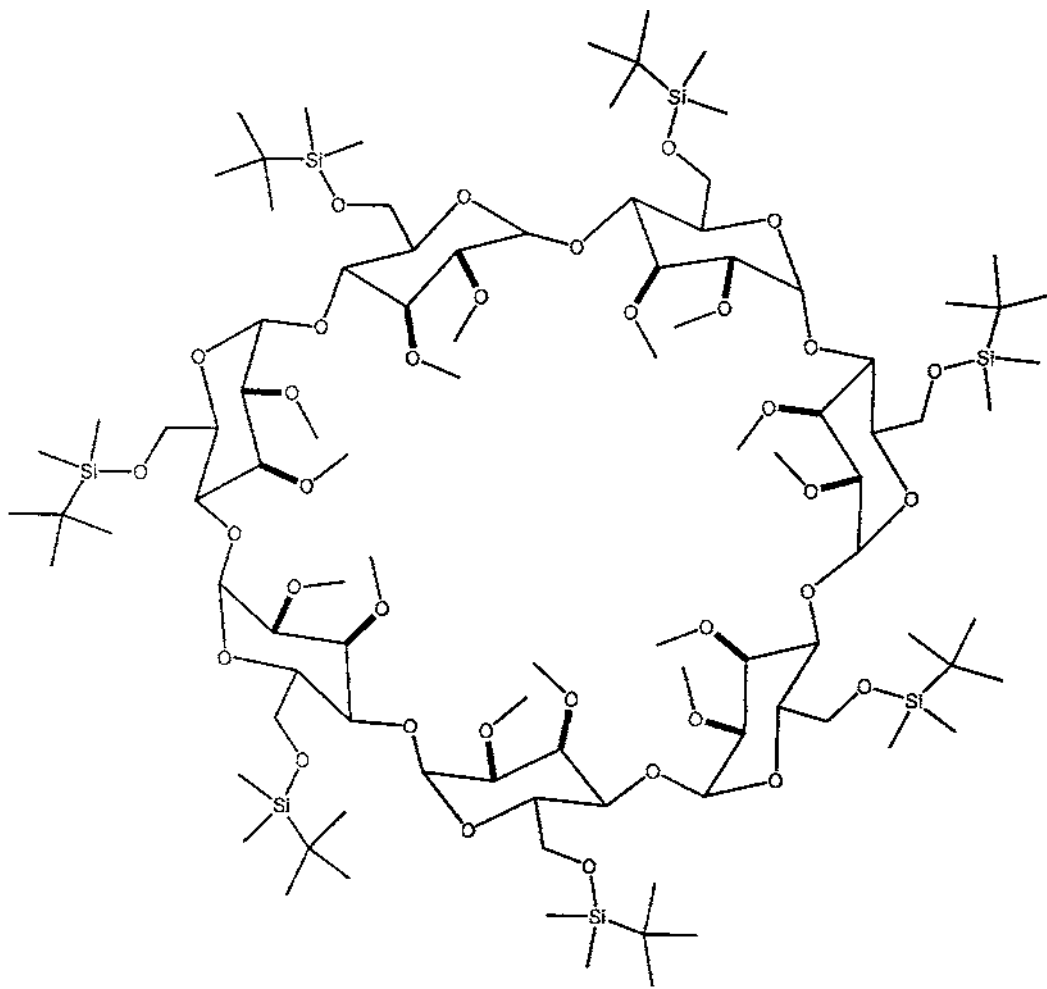


Figure 22.7
Cyclodextrin derivative used as a chiral adsorbent in gas chromatography.

this method was used to isolate a pure sample of (S)-(+)-2-methyl-4-octanol, a male-specific compound released by the sugarcane weevil, *Sphenophorus levis*.²

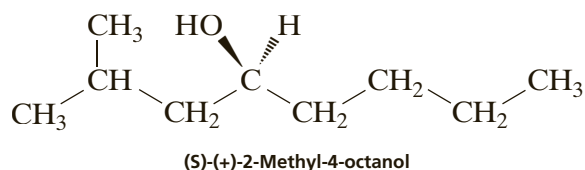

22.9 Poor Resolution and Tailing

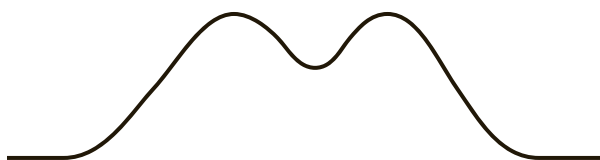


Figure 22.8
Poor resolution or peaks overlap.

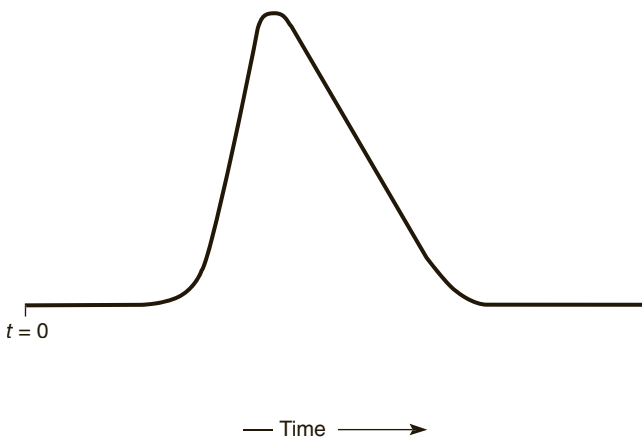


Figure 22.9
Tailing.

22.10 Qualitative Analysis

The peaks in Figure 22.5 are well **resolved**. That is, the peaks are separated from one another, and between each pair of adjacent peaks the tracing returns to the **baseline**.

In Figure 22.8, the peaks overlap, and the resolution is not good. Poor resolution is often caused by using too much sample, too high a column temperature, too short a column, a liquid phase that does not discriminate well between the two components, a column with too large a diameter, or, in short, almost any wrongly adjusted parameter. When peaks are poorly resolved, it is more difficult to determine the relative amount of each component. Methods for determining the relative percentages of each component are given in Section 22.12.

Another desirable feature illustrated by the chromatogram in Figure 22.5 is that each peak is symmetrical. A common example of an unsymmetrical peak is one in which **tailing** has occurred, as shown in Figure 22.9. Tailing usually results from injecting too much sample into the gas chromatograph. Another cause of tailing occurs with polar compounds, such as alcohols and aldehydes. These compounds may be temporarily adsorbed on column walls or areas of the support material that are not adequately coated by the liquid phase. Therefore, they do not leave in a band, and tailing results.

A disadvantage of the gas chromatograph is that it gives no information whatever about the identities of the substances it has separated. The little information it does provide is given by the retention time. It is hard to reproduce this quantity from day to day, however, and exact duplications of separations performed last month may be difficult to make this month. It is usually necessary to **calibrate** the column each time it is used. That is, you must run pure samples of all known and suspected components of a mixture individually, just before chromatographing the mixture, to obtain the retention time of each known compound. As an alternative, each suspected component can be added, one by one, to the unknown mixture while the operator looks to see which peak has its intensity increased relative to the unmodified mixture. Another solution is to *collect* the components individually as they emerge from the gas chromatograph. Each component can then be identified by other means, such as by infrared or nuclear magnetic resonance spectroscopy or by mass spectrometry.

² Zarbin, P. H. G.; Princival, J. L.; dos Santos, A. A.; de Oliveira, A. R. M. Synthesis of (S)-(+)-2-Methyl-4-octanol: Male-Specific Compound Released by Sugarcane Weevil *Sphenophorus levis*. *J. Braz. Chem. Soc.* 15 (2004): 331–334.

22.11 Collecting the Sample

For gas chromatography with a thermal conductivity detector, it is possible to collect samples that have passed through the column. One method uses a gas-collection tube (see Figure 22.10), which is included in most microscale glassware kits. A collection tube is joined to the exit port of the column by inserting the $\text{F}5/5$ inner joint into a metal adapter, which is connected to the exit port. When a sample is eluted from the column in the vapor state, it is cooled by the connecting adapter and the gas-collection tube and condenses in the collection tube. The gas-collection tube is removed from the adapter when the recorder indicates that the desired sample has completely passed through the column. After the first sample has been collected, the process can be repeated with another gas-collection tube.

To isolate the liquid, the tapered joint of the collection tube is inserted into a 0.1-mL conical vial, which has a $\text{F}5/5$ outer joint. The assembly is placed into a test tube, as illustrated in Figure 22.11. During centrifugation, the sample is forced into the bottom of the conical vial. After the apparatus is disassembled, the liquid can be removed from the vial with a syringe for a boiling-point determination or analysis by infrared spectroscopy. If a determination of the sample weight is desired, the empty conical vial and cap should be tared and reweighed after the liquid has been collected. It is advisable to dry the gas collection tube and the conical vial in an oven before use to prevent contamination by water or other solvents used in cleaning this glassware.

Another method for collecting samples is to connect a cooled trap to the exit port of the column. A simple trap, suitable for microscale work, is illustrated in Figure 22.12. Suitable coolants include ice water, liquid nitrogen, or

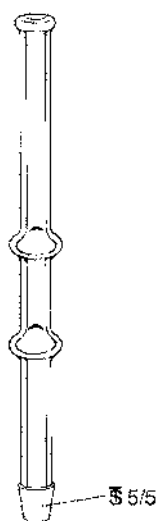


Figure 22.10
Gas-chromatography collection tube.

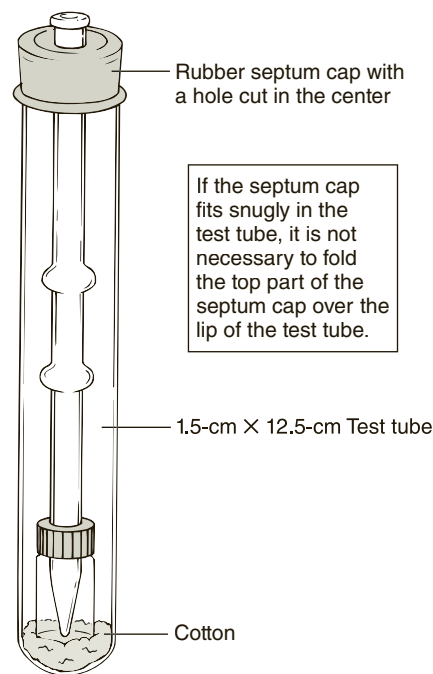


Figure 22.11
Gas-chromatography collection tube and 0.1-mL conical vial.

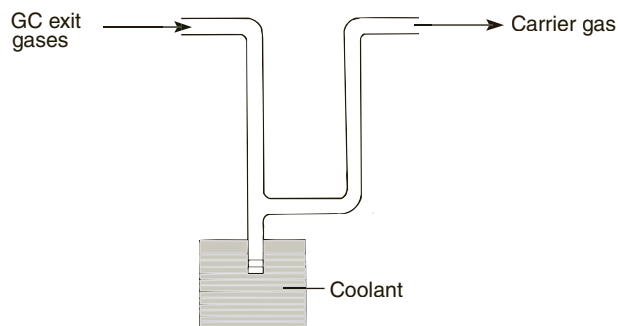


Figure 22.12
Collection trap.

dry ice–acetone. For instance, if the coolant is liquid nitrogen (bp -196°C) and the carrier gas is helium (bp -269°C), compounds boiling above the temperature of liquid nitrogen generally are condensed or trapped in the small tube at the bottom of the U-shaped tube. The small tube is scored with a file just below the point where it is connected to the larger tube, the tube is broken off, and the sample is removed for analysis. To collect each component of the mixture, you must change the trap after each sample is collected.

22.12 Quantitative Analysis

The area under a gas-chromatograph peak is proportional to the amount (moles) of compound eluted. Hence, the molar percentage composition of a mixture can be approximated by comparing relative peak areas. This method of analysis assumes that the detector is equally sensitive to all compounds eluted and that it gives a linear response with respect to amount. Nevertheless, it gives reasonably accurate results.

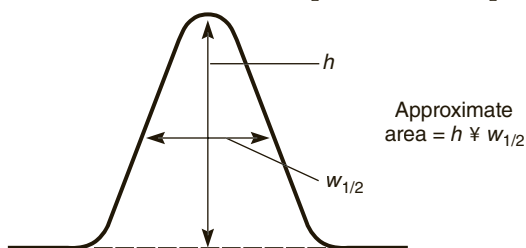


Figure 22.13
Triangulation of a peak.

The simplest method of measuring the area of a peak is by geometrical approximation, or triangulation. In this method, you multiply the height h of the peak above the baseline of the chromatogram by the width of the peak at half of its height $w_{1/2}$. This is illustrated in Figure 22.13. The baseline is approximated by drawing a line between the two side arms of the peak. This method works well only if the peak is symmetrical. If the peak has tailed or is unsymmetrical, it is best to cut out the peaks with scissors and weigh the pieces of paper on an **analytical balance**. Because the weight per area of a piece of good chart paper is reasonably constant from place to place, the ratio of the areas is the same as the ratio of the weights.

To obtain a percentage composition for the mixture, first add all the peak areas (weights). Then, to calculate the percentage of any component in the mixture, divide its individual area by the total area and multiply the result by 100. A sample calculation is illustrated in Figure 22.14. If peaks overlap (see Figure 22.8), either the gas-chromatographic conditions must be readjusted to achieve better resolution of the peaks or the peak shape must be estimated.

There are various instrumental means, which are built into recorders, of detecting the amounts of each sample automatically. One method uses a separate pen that produces a trace that integrates the area under each peak. Another method employs an electronic device that automatically prints out the area under each peak and the percentage composition of the sample.

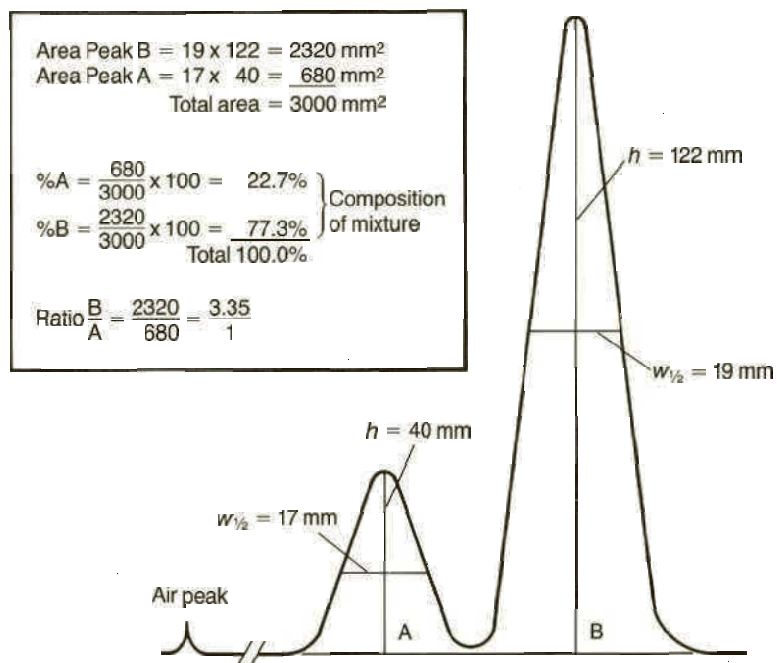


Figure 22.14
Sample percentage composition calculation.

Most modern data stations label the top of each peak with its retention time in minutes. When the trace is completed, the computer prints a table of all the peaks with their retention times, areas, and the percentage of the total area (sum of all the peaks) that each peak represents. Some caution should be used with these results because the computer often does not include smaller peaks and occasionally does not resolve narrow peaks that are so close together that they overlap. If the trace has several peaks and you would like the ratio of only two of them, you will have to determine their percentages yourself using only their two areas or instruct the instrument to integrate only these two peaks.

For the experiments in this textbook, we have assumed that the detector is equally sensitive to all compounds eluted. Compounds with different functional groups or with widely varying molecular weights, however, produce different responses with both TCD and FID gas chromatographs. With a TCD, the responses are different because not all compounds have the same thermal conductivity. Different compounds analyzed with an FID gas chromatograph also give different responses because the detector response varies with the type of ions produced. For both types of detectors, it is possible to calculate a **response factor** for each compound in a mixture. Response factors are usually determined by making up an equimolar mixture of two compounds, one of which is considered to be the reference. The mixture is separated on a gas chromatograph, and the relative percentages are calculated using one of the methods described previously. From these percentages you can determine a response factor for the compound being compared to the reference. If you do this for all the components in a mixture, you can then use these correction factors to make more accurate calculations of the relative percentages for the compounds in the mixture.

To illustrate how response factors are determined, consider the following example. An equimolar mixture of benzene, hexane, and ethyl acetate is

prepared and analyzed using a flame-ionization gas chromatograph. The peak areas obtained are

Hexane 831158
Ethyl acetate 1449695
Benzene 966463

In most cases, benzene is taken as the standard, and its response factor is defined to be equal to 1.00. Calculation of the response factors for the other components of the test mixture proceeds as follows:

Hexane $831158/966463 = 0.86$
Ethyl acetate $1449695/966463 = 1.50$
Benzene $966463/966463 = 1.00$ (by definition)

Notice that the response factors calculated in this example are molar response factors. It is necessary to correct these values by the relative molecular weights of each substance to obtain weight response factors.

When you use a flame-ionization gas chromatograph for quantitative analysis, it is first necessary to determine the response factors for each component of the mixture being analyzed, as just shown. For a quantitative analysis, it is likely that you will have to convert molar response factors into weight response factors. Next, the chromatography experiment using the unknown samples is performed. The observed peak areas for each component are corrected using the response factors in order to arrive at the correct weight percentage of each component in the sample. The application of response factors to correct the original results of a quantitative analysis will be illustrated in the following section.

22.13 Treatment of Data: Chromatograms Produced by Modern Data Stations

A. Gas Chromatograms and Data Tables

Most modern gas chromatography instruments are equipped with computer-based data stations. Interfacing the instrument with a computer allows the operator to display and manipulate the results in whatever manner might be desired. The operator thus can view the output in a convenient form. The computer can display the actual gas chromatogram and display the integration results. It can even display the result of two experiments simultaneously, making a comparison of parallel experiments convenient.

Figure 22.15 shows a gas chromatogram of a mixture of hexane, ethyl acetate, and benzene. The peaks corresponding to each substance can be seen; the peaks are labeled with their respective retention times. We can also see that there is a very small amount of an unspecified impurity, with a retention time of about 3.4 minutes.

<i>Retention Time (minutes)</i>	
Hexane	2.959
Ethyl acetate	3.160
Benzene	3.960

Figure 22.16 shows part of the printed output that accompanies the gas chromatogram. It is this information that is used in the quantitative analysis of the mixture. According to the printout, the first peak has a retention time of 2.954 minutes (the difference between the retention times that appear as labels on the graph and those that appear in the data table are not significant). The computer has also determined the area under this peak (422,373 counts). Finally, the computer has calculated the percentage of the first substance (hexane) by determining the total area

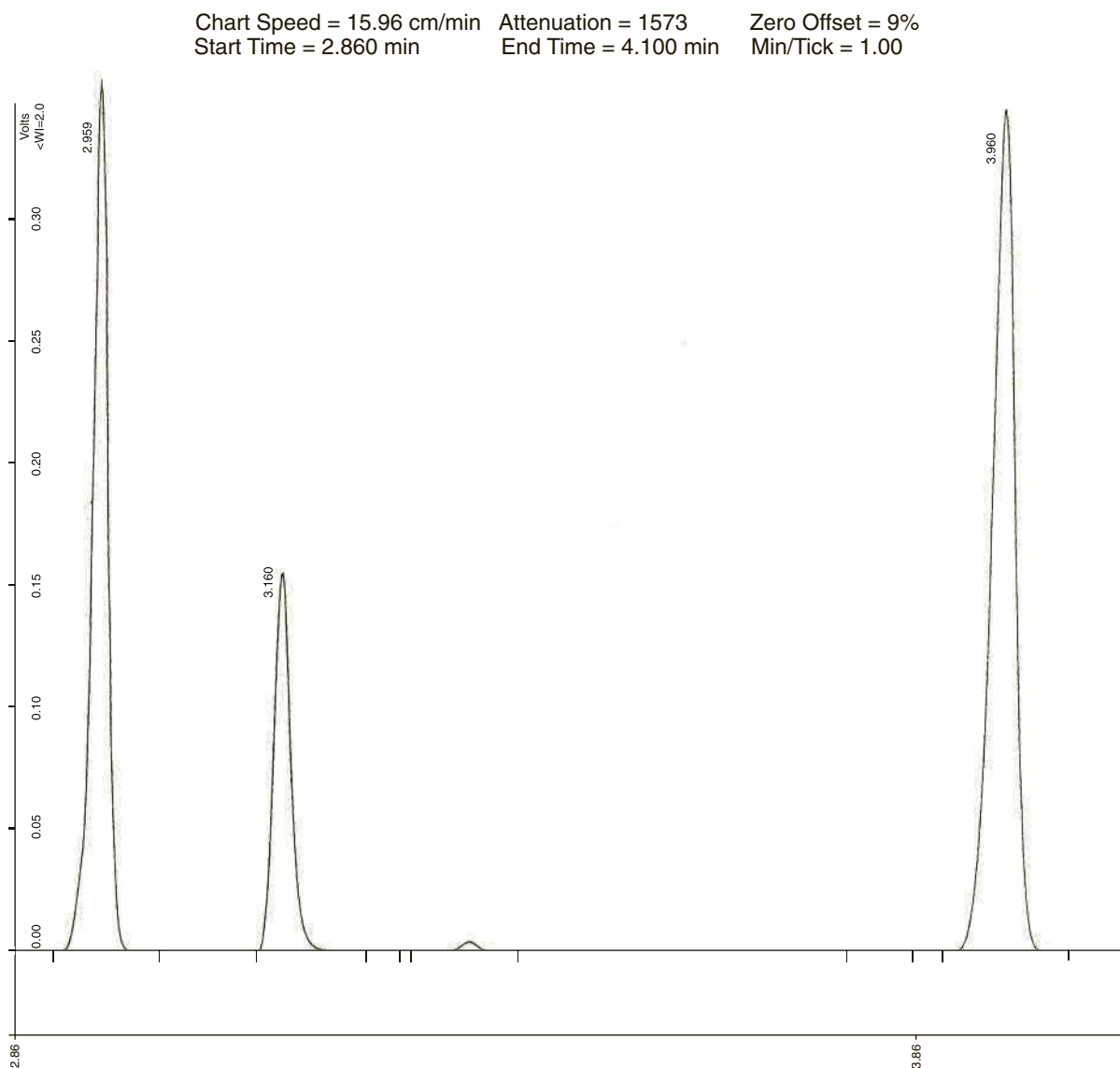


Figure 22.15
A sample gas chromatogram obtained from a data station.

of all the peaks in the chromatogram (1,227,054 counts) and dividing that into the area for the hexane peak. The result is displayed as 34.4217%. In a similar manner, the data table shows the retention times and peak areas for the other two peaks in the sample, along with a determination of the percentage of each substance in the mixture.

B. Application of Response Factors

If the detector responded with equal sensitivity to each of the components of the mixture, the data table shown in Figure 22.16 would contain the complete quantitative analysis of the sample. Unfortunately, as we have seen (Section 22.12), gas chromatography detectors respond more sensitively to some substances than they do to others. To correct for this discrepancy, it is necessary to apply corrections that are based on the **response factors** for each component of the mixture.

The method for determining the response factors was introduced in Section 22.12. In this section, we will see how this information is applied in order to obtain a correct

Run Mode : Analysis
 Peak Measurement: Peak Area
 Calculation Type: Percent

Peak No.	Peak Name	Result ()	Ret. Time (min)	Time Offset (min)	Area (counts)	Sep. Code	Width 1/2 (sec)	Status Codes
1		34.4217	2.954	0.000	422373	BB	1.0	
2		16.6599	3.155	0.000	204426	BB	1.2	
3		48.9184	3.954	0.000	600255	BB	1.6	
Totals:		100.0000		0.000	1227054			

Total Unidentified Counts : 1227054 counts

Detected Peaks: 8 Rejected Peaks: 5 Identified Peaks: 0

Multiplier: 1 Divisor: 1 Unidentified Peak Factor: 0

Baseline Offset: 1 microVolts

Noise (used): 28 microVolts — monitored before this run

Manual injection

Figure 22.16

A data table to accompany the gas chromatogram shown in Figure 22.15.

analysis. This example should serve to demonstrate the procedure for correcting raw gas chromatography results when response factors are known. According to the data table, the reported peak area for the first (hexane) peak is 422,373 counts. The response factor for hexane was previously determined to be 0.86. The area of the hexane peak is thus corrected as follows:

$$422,373/0.86 = 491,000$$

Notice that the calculated result has been adjusted to reflect a reasonable number of significant figures.

The areas for the other peaks in the gas chromatogram are corrected in a similar manner:

Hexane	422,373/0.86 =	491,000
Ethyl acetate	204,426/1.50 =	136,000
Benzene	600,255/1.00 =	600,000
Total peak area		1,227,000

Using these corrected areas, the true percentages of each component can be easily determined:

		Composition
Hexane	491,000/1,227,000	40.0%
Ethyl acetate	136,000/1,227,000	11.1%
Benzene	600,000/1,227,000	48.9%
Total		100.0%

C. Determination of Relative Percentages of Components in a Complex Mixture

In some circumstances, one may wish to determine the relative percentages of two components when the mixture being analyzed may be more complex and may contain more than two components. Examples of this situation might include the analysis of a reaction product where the laboratory worker might be interested in the relative percentages of two isomeric products when the sample might also contain peaks arising from the solvent, unreacted starting material, or some other product or impurity.

The example provided in Figures 22.15 and 22.16 can be used to illustrate the method of determining the relative percentages of some, but not all, of the components in the sample. Assume we are interested in the relative percentages of hexane and ethyl acetate in the sample but not in the percentage of benzene, which may be a solvent or an impurity. We know from the previous discussion that the *corrected* relative areas of the two peaks of interest are as follows:

	<i>Relative Area</i>
Hexane	491,000
Ethyl acetate	136,000
Total	627,000

We can determine the relative percentages of the two components simply by dividing the area of each peak by the total area of the two peaks:

		<i>Percentage</i>
Hexane	491,000/627,000	78.3%
Ethyl acetate	136,000/627,000	21.7%
Total		100.0%

22.14 Gas Chromatography—Mass Spectrometry (GC-MS)

A variation on gas chromatography is **gas chromatography-mass spectrometry**, also known as **GC-MS**. In this technique, a gas chromatograph is coupled to a mass spectrometer (see Technique 28). In effect, the mass spectrometer acts as a detector. The gas stream emerging from the gas chromatograph is admitted through a valve into a tube, where it passes over the sample inlet system of the mass spectrometer. Some of the gas stream is thus admitted into the ionization chamber of the mass spectrometer.

The molecules in the gas stream are converted into ions in the ionization chamber, and thus the gas chromatogram is actually a plot of time versus **ion current**, a measure of the number of ions produced. At the same time that the molecules are converted into ions, they are also accelerated and passed through the **mass analyzer** of the instrument. The instrument, therefore, determines the mass spectrum of each fraction eluting from the gas chromatography column.

A drawback of this method involves the need for rapid scanning by the mass spectrometer. The instrument must determine the mass spectrum of each component in the mixture before the next component exits from the column so that the spectrum of one substance is not contaminated by the spectrum of the next fraction.

Because high-efficiency capillary columns are used in the gas chromatograph, in most cases compounds are completely separated before the gas stream is analyzed. The typical GC-MS instrument has the capability of obtaining at least one scan per second in the range of 10–300 amu. Even more scans are possible if a narrow range of masses is analyzed. Using capillary columns, however, requires the user to take particular care to ensure that the sample does not contain any particles that might obstruct

the flow of gases through the column. For this reason, the sample is carefully filtered through a very fine filter before the sample is injected into the chromatograph.

With a GC–MS system, a mixture can be analyzed and results obtained that resemble very closely those shown in Figures 22.15 and 22.16. A library search on each component of the mixture can also be conducted. The data stations of most instruments contain a library of standard mass spectra in their computer memory. If the components are known compounds, they can be identified tentatively by a comparison of their mass spectrum with the spectra of compounds found in the computer library. In this way, a “hit list” can be generated that reports on the probability that the compound in the library matches the known substance. A typical printout from a GC–MS instrument will list probable compounds that fit the mass spectrum of the component, the names of the compounds, their CAS Nos. (see Technique 29, Section 29.11), and a “quality” or “confidence” number. This last number provides an estimate of how closely the mass spectrum of the component matches the mass spectrum of the substance in the computer library.

A variation on the GC–MS technique includes coupling a Fourier transform infrared spectrometer (FT–IR) to a gas chromatograph. The substances that elute from the gas chromatograph are detected by determining their infrared spectra rather than their mass spectra. A new technique that also resembles GC–MS is **high-performance liquid chromatography–mass spectrometry (HPLC–MS)**. An HPLC instrument is coupled through a special interface to a mass spectrometer. The substances that elute from the HPLC column are detected by the mass spectrometer, and their mass spectra can be displayed, analyzed, and compared with standard spectra found in the computer library built into the instrument.

PROBLEMS

- A sample consisting of 1-bromopropane and 1-chloropropane is injected into a gas chromatograph equipped with a nonpolar column. Which compound has the shorter retention time? Explain your answer.
 - If the same sample were run several days later with the conditions as nearly the same as possible, would you expect the retention times to be identical to those obtained the first time? Explain.
- Using triangulation, calculate the percentage of each component in a mixture composed of two substances, A and B. The chromatogram is shown in Figure 22.17.

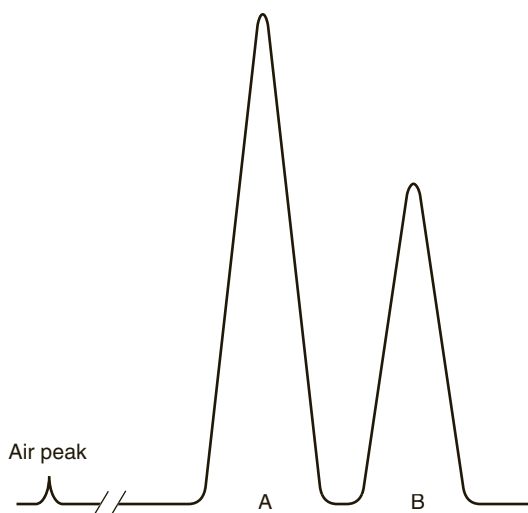


Figure 22.17
A chromatogram for problem 2.

3. Make a photocopy of the chromatogram in Figure 22.17. Cut out the peaks and weigh them on an analytical balance. Use the weights to calculate the percentage of each component in the mixture. Compare your answer to what you calculated in problem 2.
4. What would happen to the retention time of a compound if the following changes were made?
 - a. Decrease the flow rate of the carrier gas
 - b. Increase the temperature of the column
 - c. Increase the length of the column

Polarimetry

23.1 Nature of Polarized Light

Light has a dual nature because it shows properties of both waves and particles. The wave nature of light can be demonstrated by two experiments: polarization and interference. Of the two, polarization is the more interesting to organic chemists because they can take advantage of polarization experiments to learn something about the structure of an unknown molecule.

Ordinary white light consists of wave motion in which the waves have a variety of wavelengths and vibrate in all possible planes perpendicular to the direction of propagation. Light can be made to be **monochromatic** (of one wavelength or color) by using filters or special light sources. Frequently, a sodium lamp (sodium D line = 5893 Å) is used. Although the light from this lamp consists of waves of only one wavelength, the individual light waves still vibrate in all possible planes perpendicular to the beam. If we imagine that the beam of light is aimed directly at the viewer, ordinary light can be represented by showing the edges of the planes oriented randomly around the path of the beam, as on the left side of Figure 23.1.

A Nicol prism, which consists of a specially prepared crystal of Iceland spar (or calcite), has the property of serving as a screen that can restrict the passage of light waves. Waves that are vibrating in one plane are transmitted; those in all other planes are rejected (either refracted in another direction or absorbed). The light that passes through the prism is called **plane-polarized light**, and it consists of waves that vibrate in only one plane. A beam of plane-polarized light aimed directly at the viewer can be represented by showing the edges of the plane oriented in one particular direction, as on the right side of Figure 23.1.

Iceland spar has the property of **double refraction**; that is, it can split, or doubly refract, an entering beam of ordinary light into two separate emerging beams of light. Each of the two emerging beams (labeled A and B in Figure 23.2) has only a single plane of vibration, and the plane of vibration in beam A is perpendicular to the plane of vibration in beam B. In other words, the crystal has separated the incident beam of

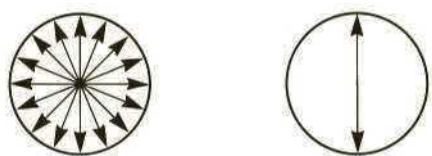


Figure 23.1
Ordinary versus plane-polarized light.

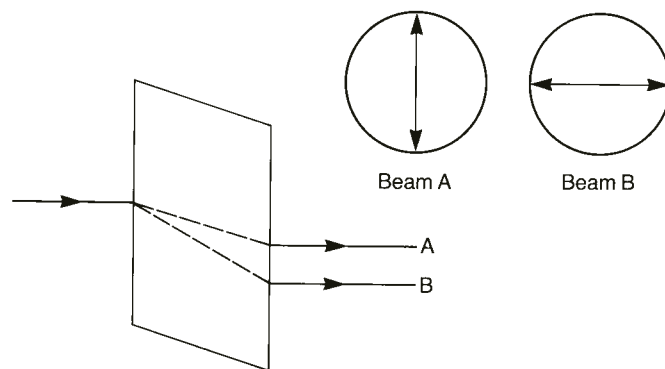


Figure 23.2
Double refraction.

ordinary light into two beams of plane-polarized light, with the plane of polarization of beam A perpendicular to the plane of beam B.

To generate a single beam of plane-polarized light, one can take advantage of the double-refracting property of Iceland spar. A Nicol prism, invented by the Scottish physicist William Nicol, consists of two crystals of Iceland spar cut to specified angles and cemented by Canada balsam. This prism transmits one of the two beams of plane-polarized light while reflecting the other at a sharp angle so that it does not interfere with the transmitted beam. Plane-polarized light can also be generated by a Polaroid filter, a device invented by E. H. Land, an American physicist. Polaroid filters consist of certain types of crystals embedded in transparent plastic and capable of producing plane-polarized light.

After passing through a first Nicol prism, plane-polarized light can pass through a second Nicol prism, but only if the second prism has its axis oriented so that it is *parallel* to the incident light's plane of polarization. Plane-polarized light is *absorbed* by a Nicol prism that is oriented so that its axis is *perpendicular* to the incident light's plane of polarization. These situations can be illustrated by the picket-fence analogy, as shown in Figure 23.3. Plane-polarized light can pass through a fence whose slats are oriented in the proper direction but is blocked out by a fence whose slats are oriented perpendicularly.

An **optically active substance** is one that interacts with polarized light to rotate the plane of polarization through some angle α . Figure 23.4 illustrates this phenomenon.

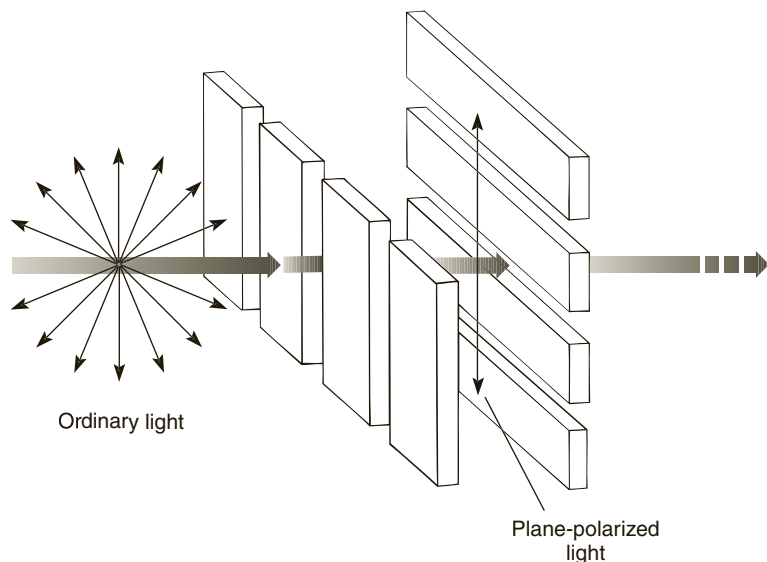


Figure 23.3
The picket-fence analogy.

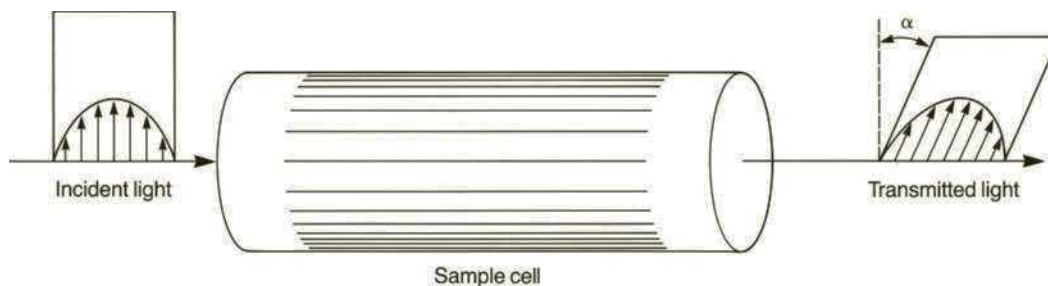


Figure 23.4
Optical activity.

23.2 The Polarimeter An instrument called a **polarimeter** is used to measure the extent to which a substance interacts with polarized light. A schematic diagram of a polarimeter is shown in Figure 23.5. The light from the source lamp is polarized by being passed through a fixed Nicol prism, called a **polarizer**. This light passes through the sample, with which it may or may not interact to have its plane of polarization rotated in one direction or the other. A second, rotatable Nicol prism, called the **analyzer**, is adjusted to allow the maximum amount of light to pass through. The number of degrees and the direction of rotation required for this adjustment are measured to give the observed rotation α .

So that data determined by several persons under different conditions can be compared, a standardized means of presenting optical rotation data is necessary. The most common way of presenting such data is by recording the **specific rotation** $[\alpha]_{\lambda}^t$, which has been corrected for differences in concentration, cell path length, temperature, solvent, and wavelength of the light source. The equation defining the specific rotation of a compound in solution is

$$[\alpha]_{\lambda}^t = \frac{\alpha}{cl}$$

where α = observed rotation in degrees, c = concentration in grams per milliliter of solution, l = length of sample tube in decimeters, λ = wavelength of light (usually indicated as "D," for the sodium D line), and t = temperature in degrees Celsius. For pure liquids, the density d of the liquid in grams per milliliter replaces c in the preceding formula. You may occasionally want to compare compounds of different molecular weights, so a **molecular rotation**, based on moles instead of grams, is more convenient than a specific rotation.

The molecular rotation M_{λ}^t is derived from the specific rotation $[\alpha]_{\lambda}^t$ by

$$M_{\lambda}^t = \frac{[\alpha]_{\lambda}^t \times \text{Molecular weight}}{100}$$

Usually, measurements are made at 25°C with the sodium D line as a light source; consequently, specific rotations are reported as $[\alpha]_{\text{D}}^{25}$.

Polarimeters that are now available incorporate electronics to determine the angle of rotation of chiral molecules. These instruments are essentially automatic. The only real difference between an automatic polarimeter and a manual one is that a light detector replaces the eye. No visual observation of any kind is made with an automatic instrument. A microprocessor adjusts the analyzer until the light reaching the detector is at a minimum. The angle of rotation is displayed digitally in an LCD window, including the sign of rotation. The simplest instrument is equipped with a sodium lamp that gives rotations based on the sodium D line (589 nm). More

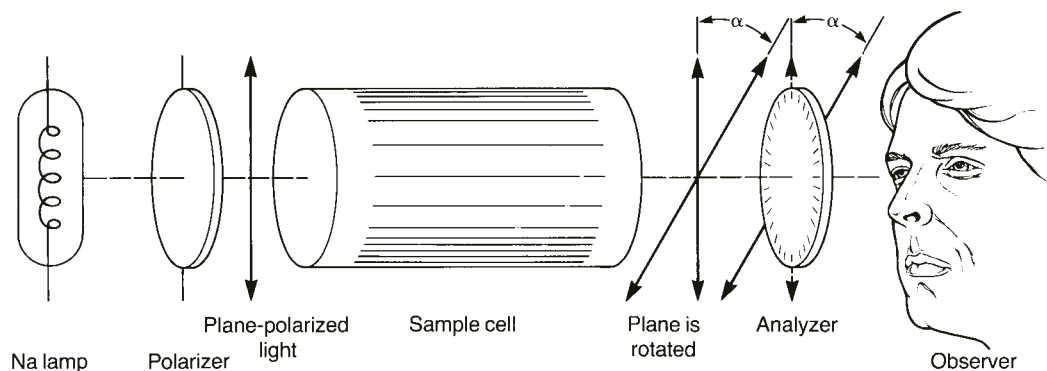


Figure 23.5
Schematic diagram of a polarimeter.

expensive instruments use a tungsten lamp and filters so that wavelengths can be varied over a range of values. Using the latter instrument, a chemist can observe rotations at different wavelengths.

23.3 Sample Preparation: The Sample Cell

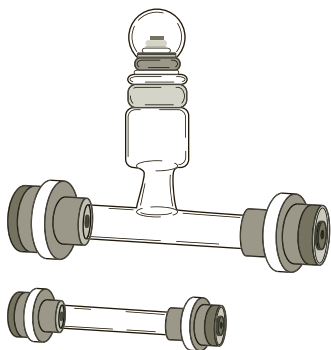


Figure 23.6
Two modern polarimetry cells (Rudolph Research).

It is important for the solution whose optical rotation is to be determined to contain no suspended particles of dust, dirt, or undissolved material that might disperse the incident polarized light. Therefore, you must clean the sample cell carefully, and your sample must be free of suspended particles. You must also prevent the presence of any air bubbles in the bore when you fill the cell. Most cells have a stem in the center or an area at one end of the cell where the diameter of the tube is increased. These features are designed to help you catch any bubbles in an area that is above the path that the light takes through the main bore.

Two modern **polarimetry cells** are shown in Figure 23.6. In the first case, the cell is filled until the liquid completely fills the bore and a small portion of the center stem. Then, if one gently rocks the cell back and forth along its axis, bubbles will rise and collect in the stem where they are above the light path. A stopper is placed in the stem when you are finished. In the second case, the cell is filled vertically, and the end is screwed on. Bubbles are trapped at the raised end when the cell is turned horizontally.

Sample cells are available in various lengths, with 0.5 dm and 1.0 dm being the most common. A typical 0.5-dm cell holds about 3–5 mL of solution, but many companies sell **microcells** that have a very narrow-diameter bore and require much less solution. Polarimeter cells are quite expensive because the windows must be made out of quartz rather than ordinary glass. Be sure to handle them carefully and to avoid getting fingerprints on the end windows because this will also disperse the polarized light.

With liquid samples, it is often possible to use the **neat** (undiluted) liquid as your sample. In this case, the concentration of the sample is just the density of the liquid (g/mL). If you have a solid sample or if you have too little of a liquid to fill the cell, you will have to either dissolve or dilute the sample with a solvent. In this case, you must weigh (grams) the amount of material you use and divide by the total volume (mL) to obtain the concentration in g/mL. Water, methanol, and ethanol are the best solvents to use because they are unlikely to attack the cell you are using. Many cells have rubber parts or use a cement to attach the windows to the ends of the bore. Rubber and cements will often dissolve in stronger solvents such as acetone or methylene chloride, thereby damaging the cell. Check with your instructor before using any solvent stronger than water, methanol, or ethanol. These are also the preferred solvents to use for cleaning the cells.

23.4 Operation of the Polarimeter

A. The Zeiss Polarimeter, a Classic Instrument

The procedures given here are for the operation of the Zeiss polarimeter (Figure 23.7), a classic analog instrument with a circular scale and a sodium lamp. Many other older models of polarimeter are operated in a similar fashion.

Before taking any measurements, turn on the sodium lamp and wait 5–10 minutes for the lamp to warm up and stabilize. After the warm-up period is complete, you should make an initial check of the instrument by taking a zero reading with a sample cell filled only with solvent. If the zero reading does not correspond with the zero-degree (0°) calibration mark, then the difference in readings must be used to correct all subsequent readings.

To take the zero measurement, place the polarimeter cell with the sample in the sloped cradle or rack inside the instrument. If you are using a cell with an enlarged end, that end must be placed at the high end of the cradle, making sure that no

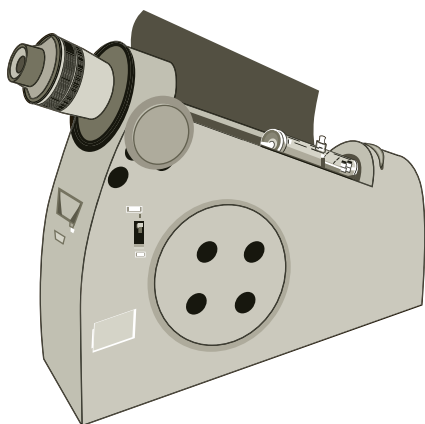



Figure 23.7
The Zeiss polarimeter.

bubbles are in the bore of the cell. After closing the cover and while watching through the eyepiece, turn the analyzer knob or ring until the proper angle of the analyzer is reached (the angle that allows no light to pass through the instrument). Most analog instruments, including the Zeiss polarimeter, are of the split-field type. When you look upward through the eyepiece, you see a circle split into three sectors (Figure 23.8), with the center sector either lighter or darker than those on either side. The analyzer prism is rotated until all of the sectors are matched in intensity, usually the darker color (see Figure 23.8). This is called the **null** reading.

When you look downward in the eyepiece, you see the value of the angle through which the plane of the polarized light has been rotated (if any) indicated on a vernier degree scale (Figure 23.9). Some polarimeters, such as the original Rudolph polarimeter, have instead a large circular scale, like a halo, attached directly to the knob you turn.


After determining the zero setting on the blank solution, place the polarimeter cell containing your sample into the polarimeter and measure the observed angle of rotation in the same way as described for the zero measurement. Be sure to record not only the numerical value of the reading but also the direction of rotation. Also record the solvent, temperature, and concentration, as these are also critical to the measurement. Rotations clockwise are due to dextrorotatory substances and are indicated by the “+” sign. Rotations counterclockwise are due to levorotatory substances and are indicated by the “-” sign. You should take several readings, including readings for which the value was approached from both sides. In other words, where the actual reading might be $+75^\circ$, first approach this reading upward from somewhere between 0° and 75° ; on the next measurement, approach the null from an angle greater than 75° . Duplicating readings, approaching the observed rotation from both sides, and averaging the readings reduce the error.

If you are not sure whether you have a dextrorotatory or a levorotatory substance, you can make this determination by halving the concentration of your compound, reducing the length of the cell by half, or reducing the intensity of the light. The confusion between dextrorotatory and levorotatory arises because you are reading a circular scale. The null reading can be approached from either direction (clockwise or counterclockwise), starting from zero (see Figure 23.10). For instance, is your null at $+120^\circ$, or is it at -240° ? Both readings are at the same point on the scale. Figure 23.10 shows that by reducing the concentration, the cell length, or the light intensity in half (any one of these), the reading will change, and it will move in a different direction for levorotatory substances than for dextrorotatory substances. The direction of rotation is most often determined by making measurements at different dilutions.

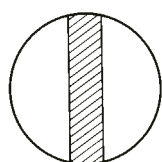


Adjustment
incorrect

“Null”

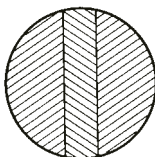

Adjustment
correct

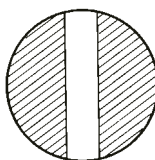

Adjustment
incorrect

Figure 23.8
Image field sectors
in the polarimeter.

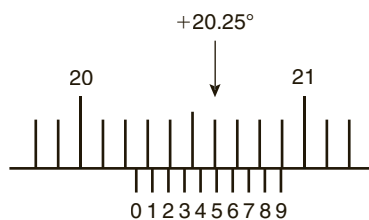


Figure 23.9
The vernier degree scale seen in the lower field
of the Zeiss polarimeter.

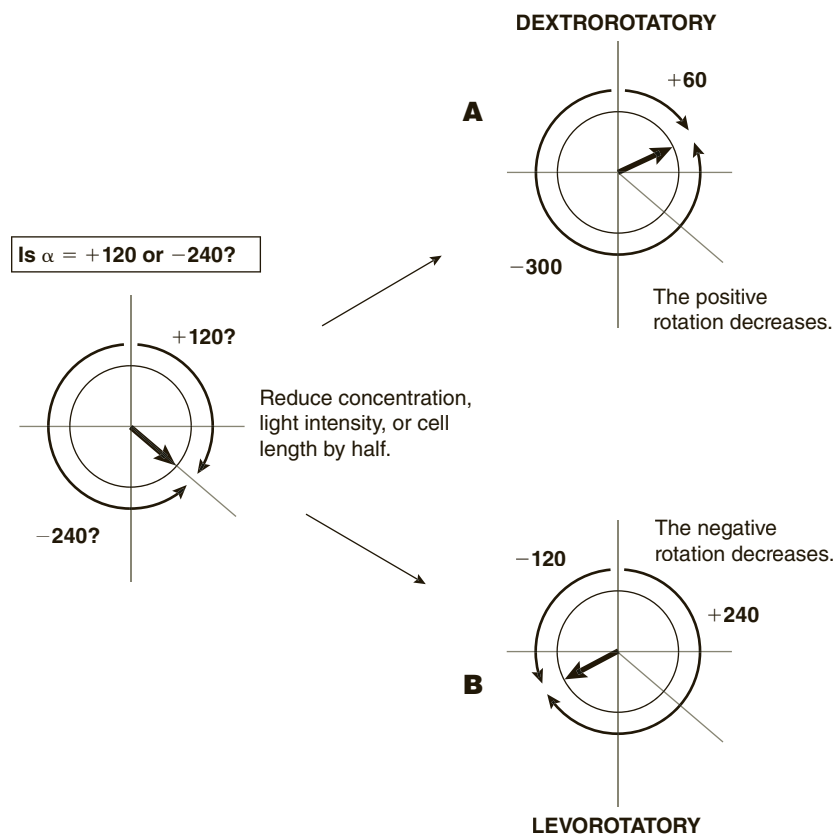


Figure 23.10

How to determine the direction of rotation. This diagram shows the effect on observed rotation if you reduce by half the concentration of the compound, the light intensity, or the length of the cell. By this method, it is easy to determine if the compound is dextrorotatory (A) or levorotatory (B).

Once you have determined the value and direction of the observed rotation α , you must correct it by the zero value and then use the formulas in Section 23.2 to convert it to the specific rotation $[\alpha]_D$. The specific rotation is always reported as a function of temperature, indicating the wavelength by "D" if a sodium lamp was used, and the solvent and concentration used are reported. For example:

$$[\alpha]_D = +43.8 \text{ (} c = 7.5 \text{ g/100 mL, in absolute ethanol)}$$

B. The Modern Digital Polarimeter

A modern digital polarimeter, such as the one shown in Figure 23.11, is much easier to operate than the older analog instruments. The modern instrument will store the zero reading for you, subtract it from every subsequent reading automatically, determine the direction of rotation, and calculate the specific rotation from the reading obtained on your sample. When finished, it can print everything on a sheet of paper for you to take with you. In a typical instrument, you first determine the zero reading and then store it in electronic memory. Once the zero reading is determined, you place your sample in the instrument. The instrument automatically finds the null angle and the direction of rotation and displays it on an LED readout. The instrument approaches the null several times to be sure of its reading and determines the direction of rotation

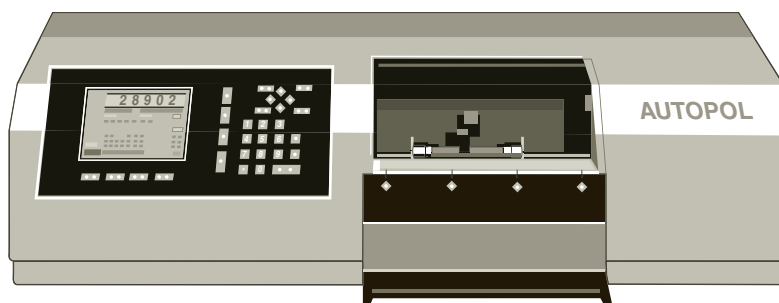


Figure 23.11
The Autopol IV (Rudolph Research), a modern digital polarimeter.

by reducing the intensity of the light. It can do this several ways. One common method is to attenuate (reduce) the incident light intensity of the beam of polarized light and see what effect this has on the angle of rotation. Even a digital polarimeter, however, cannot extract a reading from a poor sample, such as one that is cloudy, has a bubble, or has suspended solid material. A good sample is still your responsibility.

23.5 Optical Purity

When you prepare a sample of an enantiomer by a resolution method, the sample is not always 100% of a single enantiomer. It frequently is contaminated by residual amounts of the opposite stereoisomer. If you know the amount of each enantiomer in a mixture, you can calculate the **optical purity**. Some chemists prefer to use the term **enantiomeric excess (ee)** rather than optical purity. The two terms can be used interchangeably. The percentage enantiomeric excess or optical purity is calculated as follows:

$$\% \text{ Optical purity} = \frac{\text{moles one enantiomer} - \text{moles of other enantiomer}}{\text{total moles of both enantiomers}} \times 100$$

$$\% \text{ Optical purity} = \% \text{ enantiomeric excess (ee)}$$

Often, it is difficult to apply the previous equation because you do not know the exact amount of each enantiomer present in a mixture. It is far easier to calculate the optical purity (ee) by using the observed specific rotation of the mixture and dividing it by the specific rotation of the pure enantiomer. Values for the pure enantiomers can sometimes be found in literature sources.

$$\% \text{ Optical purity} = \% \text{ enantiomeric excess} = \frac{\text{observed specific rotation}}{\text{specific rotation of pure enantiomer}} \times 100$$

This latter equation holds true only for mixtures of two chiral molecules that are mirror images of each other (enantiomers). If some other chiral substance is present in the mixture as an impurity, then the actual optical purity will deviate from the value calculated.

In a racemic (\pm) mixture, there is no excess enantiomer, and the optical purity (enantiomeric excess) is zero; in a completely resolved material, the optical purity (enantiomeric excess) is 100%. A compound that is $x\%$ optically pure contains $x\%$ of one enantiomer and $(100 - x)\%$ of a racemic mixture.

Once the optical purity (enantiomeric excess) is known, the relative percentages of each of the enantiomers can be calculated easily. If the predominant form in the impure, optically active mixture is assumed to be the (+) enantiomer, the percentage of the (+) enantiomer is

$$\left[x + \left(\frac{100 - x}{2} \right) \right] \%$$

and the percentage of the (–) enantiomer is $[(100 - x)/2]\%$. The relative percentages of (+) and (–) forms in a partially resolved mixture of enantiomers can be calculated as shown next. Consider a partially resolved mixture of camphor enantiomers. The specific rotation for pure (+)-camphor is $+43.8^\circ$ in absolute ethanol, but the mixture shows a specific rotation of $+26.3^\circ$.

$$\text{Optical purity} = \frac{+26.3^\circ}{+43.8^\circ} \times 100 = 60\% \text{ optically pure}$$

$$\% (+) \text{ enantiomer} = 60 + \left(\frac{100 - 60}{2} \right) = 80\%$$

$$\% (-) \text{ enantiomer} = \left(\frac{100 - 60}{2} \right) = 20\%$$

Notice that the difference between these two calculated values equals the optical purity or enantiomeric excess.

PROBLEMS

1. Calculate the specific rotation of a substance that is dissolved in a solvent (0.4 g/mL) and that has an observed rotation of -10° as determined with a 0.5-dm cell.
2. Calculate the observed rotation for a solution of a substance (2.0 g/mL) that is 80% optically pure. A 2-dm cell is used. The specific rotation for the optically pure substance is $+20^\circ$.
3. What is the optical purity of a partially racemized product if the calculated specific rotation is -8° and the pure enantiomer has a specific rotation of -10° ? Calculate the percentage of each of the enantiomers in the partially racemized product.

Refractometry

The **refractive index** is a useful physical property of liquids. Often, a liquid can be identified from a measurement of its refractive index. The refractive index can also provide a measure of the purity of the sample being examined. This is accomplished by comparing the experimentally measured refractive index with the value reported in the literature for an ultrapure sample of the compound. The closer the measured sample's value to the literature value, the purer the sample.

24.1 The Refractive Index

The refractive index has as its basis the fact that light travels at a different velocity in condensed phases (liquids, solids) than in air. The refractive index n is defined as the ratio of the velocity of light in air to the velocity of light in the medium being measured:

$$n = \frac{V_{\text{air}}}{V_{\text{liquid}}} = \frac{\sin \theta}{\sin \phi}$$

It is not difficult to measure the ratio of the velocities experimentally. It corresponds to $(\sin \theta / \sin \phi)$, where θ is the angle of incidence for a beam of light striking the surface of the medium and ϕ is the angle of refraction of the beam of light *within* the medium. This is illustrated in Figure 24.1.

The refractive index for a given medium depends on two variable factors. First, it is *temperature* dependent. The density of the medium changes with temperature; hence, the speed of light in the medium also changes. Second, the refractive index is *wavelength* dependent. Beams of light with different wavelengths are refracted to different extents in the same medium and give different refractive indices for that

medium. It is usual to report refractive indices measured at 20°C, with a sodium discharge lamp as the source of illumination. The sodium lamp gives off yellow light of 589-nm wavelength, the so-called sodium D line. Under these conditions, the refractive index is reported in the following form:

$$n_{\text{D}}^{20} = 1.4892$$

The superscript indicates the temperature, and the subscript indicates that the sodium D line was used for the measurement. If another wavelength is used for the determination, the D is replaced by the appropriate value, usually in nanometers ($1 \text{ nm} = 10^{-9} \text{ m}$).

Notice that the hypothetical value reported has four decimal places. It is easy to determine the refractive index to within several parts in 10,000. Therefore, n_{D} is a very accurate physical constant for a given substance and can be used for identification. However, it is sensitive to even small amounts of impurity in the substance measured. Unless the substance is purified *extensively*, you will not usually be able to reproduce the last two decimal

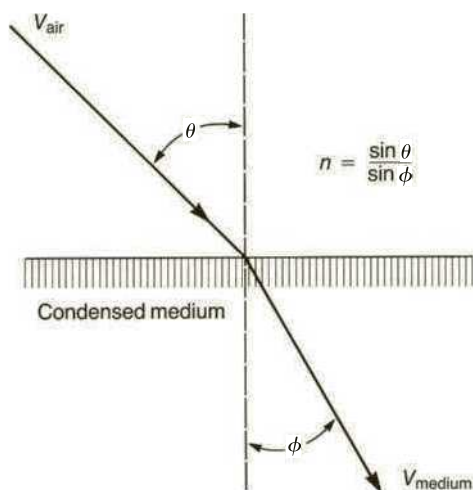


Figure 24.1
The refractive index

places given in a handbook or other literature source. Typical organic liquids have refractive index values between 1.3400 and 1.5600.

24.2 The Abbé Refractometer

The instrument used to measure the refractive index is called a **refractometer**. Although many styles of refractometer are available, by far the most common instrument is the Abbé refractometer. This style of refractometer has the following advantages:

1. White light may be used for illumination; the instrument is compensated, however, so that the index of refraction obtained is actually that for the sodium D line.
2. The prisms can be temperature controlled.
3. Only a small sample is required (a few drops of liquid using the standard method or about 5 μL using a modified technique).

A common type of Abbé refractometer is shown in Figure 24.2.

The optical arrangement of the refractometer is complex; a simplified diagram of the internal workings is given in Figure 24.3. The letters *A*, *B*, *C*, and *D* label corresponding parts in both Figures 24.2 and 24.3. A complete description of refractometer optics is too difficult to attempt here, but Figure 24.3 gives a simplified diagram of the essential operating principles.

Using the standard method, introduce the sample to be measured between the two prisms. If it is a free-flowing liquid, it may be introduced into a channel along the side of the prisms, injected from a Pasteur pipette. If it is a viscous sample, the prisms must be opened (they are hinged) by lifting the upper one; a few drops of liquid are applied to the lower prism with a Pasteur pipette or a wooden applicator. If a Pasteur pipette is used, take care not to touch the prisms because they become scratched easily. When the prisms are closed, the liquid should spread evenly to

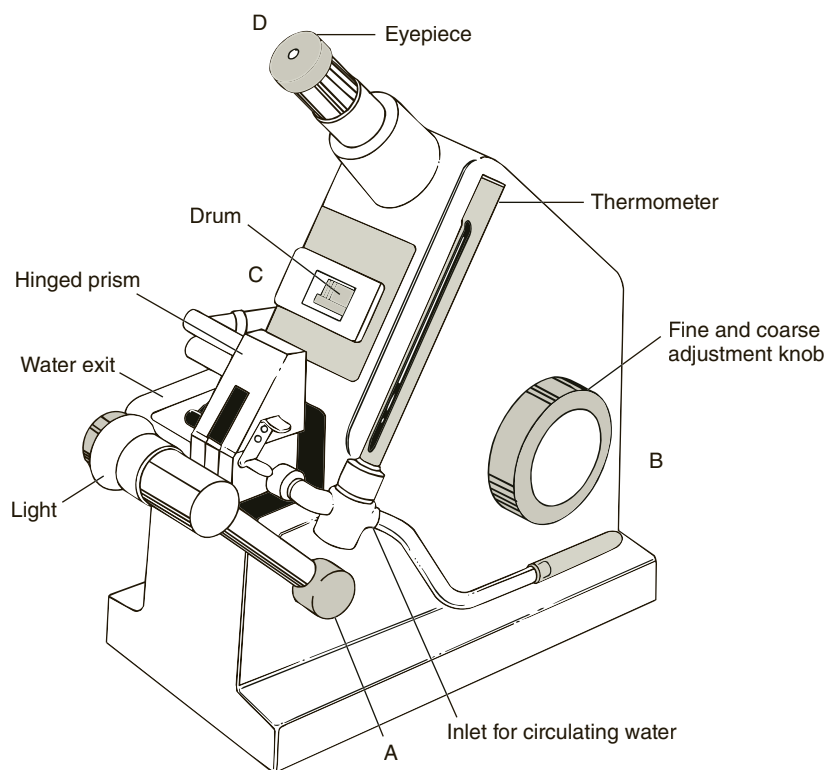


Figure 24.2
Abbé refractometer (Bausch and Lomb Abbé 3L).

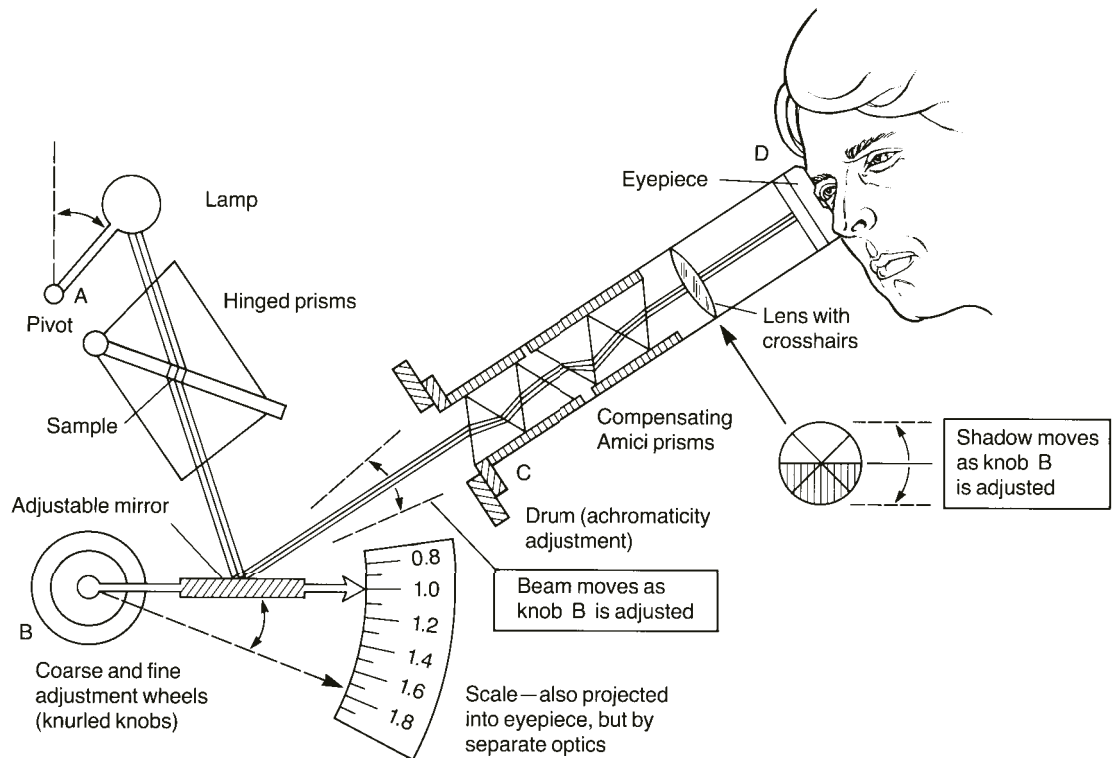


Figure 24.3
Simplified diagram of a refractometer.

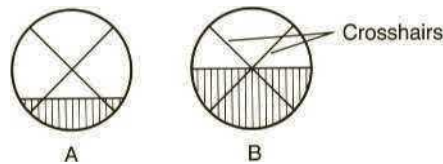


Figure 24.4
(A) Refractometer incorrectly adjusted.
(B) Correct adjustment.

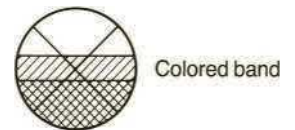


Figure 24.5
Refractometer showing chromatic aberration (color dispersion). The dispersion is incorrectly adjusted.

make a thin film. With highly volatile samples, the remaining operations must be performed rapidly. Even when the prisms are closed, evaporation of volatile liquids can readily occur.

Next, turn on the light and look into the eyepiece *D*. The hinged lamp is adjusted to give the maximum illumination to the visible field in the eyepiece. The light rotates at pivot *A*.

Rotate the coarse and fine adjustment knobs at *B* until the dividing line between the light and dark halves of the visual field coincides with the center of the crosshairs (Figure 24.4). If the crosshairs are not in sharp focus, adjust the eyepiece to focus them. If the horizontal line dividing the light and dark areas appears as a colored band, as in Figure 24.5, the refractometer shows **chromatic aberration** (color dispersion). This can be adjusted with knob *C* drum (Figure 24.3). This knurled knob rotates a series of prisms, called Amici prisms, that color-compensate the

refractometer and cancel out dispersion. Adjust the knob to give a sharp, uncolored division between the light and dark segments. When you have adjusted everything correctly (as in Figure 24.4B), read the refractive index. In the instrument described here, press a small button on the left side of the housing to make the scale visible in the eyepiece. In other refractometers, the scale is visible at all times, frequently through a separate eyepiece.

Occasionally, the refractometer will be so far out of adjustment that it may be difficult to measure the refractive index of an unknown. When this happens, it is wise to place a pure sample of known refractive index in the instrument, set the scale to the correct value of refractive index, and adjust the controls for the sharpest line possible. Once this is done, it is easier to measure an unknown sample. It is especially helpful to perform this procedure prior to measuring the refractive index of a highly volatile sample.

NOTE: There are many styles of refractometer, but most have adjustments similar to those described here.

In the procedure just described, several drops of liquid are required to obtain the refractive index. In some experiments, you may not have enough sample to use this standard method. It is possible to modify the procedure so that a reasonably accurate refractive index can be obtained on about 5 μL of liquid. Instead of placing the sample directly onto the prism, you apply the sample to a small piece of lens paper. The lens paper can be conveniently cut with a handheld paper punch,¹ and the paper disc (0.6-cm diameter) is placed in the center of the bottom prism of the refractometer. To avoid scratching the prism, use forceps or tweezers with plastic tips to handle the disc. About 5 μL of liquid is carefully placed on the lens paper using a microliter syringe. After closing the prisms, adjust the refractometer as described previously and read the refractive index. With this method, the horizontal line dividing the light and dark areas may not be as sharp as it is in the absence of the lens paper. It may also be impossible to eliminate color dispersion completely. Nonetheless, the refractive index values determined by this method are usually within 10 parts in 10,000 of the values determined by the standard procedure.

24.3 Cleaning the Refractometer

In using the refractometer, you should always remember that if the prisms are scratched, the instrument will be ruined.

NOTE: Do not touch the prisms with any hard object.

This admonition includes Pasteur pipettes and glass rods.

When measurements are completed, the prisms should be cleaned with ethanol or petroleum ether. Moisten *soft* tissues with the solvent and wipe the prisms *gently*. When the solvent has evaporated from the prism surfaces, the prisms should be locked together. The refractometer should be left with the prisms closed to avoid collection of dust in the space between them. The instrument should also be turned off when it is no longer in use.

24.4 The Digital Refractometer

Today, there are modern digital refractometers available that determine the refractive index of a liquid electronically (Figure 24.6). Once the instrument has been calibrated, it is only necessary to place a drop of your liquid between the prisms (see

¹In order to cut the lens paper more easily, place several sheets between two pieces of heavier paper, such as that used for file folders.

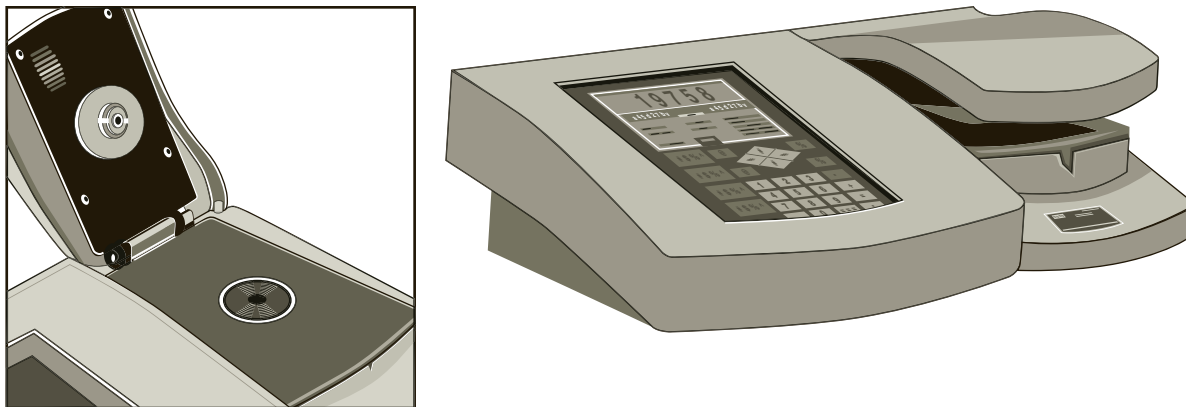


Figure 24.6

The Rudolph J-series, a modern digital refractometer. To make a measurement, place the sample on the lower prism (see the inset) and close the lid.

the inset in Figure 24.6), close the lid, and read the display. The instrument can make temperature corrections and store the values of your readings in its microprocessor memory. Once again, these instruments must be treated with respect, taking care not to scratch the prisms and to clean them after use.

24.5 Temperature Corrections

Most refractometers are designed so that circulating water at a constant temperature can maintain the prisms at 20°C. If this temperature-control system is not used or if the water is not at 20°C, a temperature correction must be made. Although the magnitude of the temperature correction may vary from one class of compound to another, a value of 0.00045 per degree Celsius is a useful approximation for most substances. The index of refraction of a substance *decreases* with *increasing* temperature. Therefore, add the correction to the observed n_D value for temperatures higher than 20°C and subtract it for temperatures lower than 20°C. For example, the reported n_D value for nitrobenzene is 1.5529. One would observe a value at 25°C of 1.5506. The temperature correction would be made as follows:

$$n_D^{20} = 1.5506 + 5(0.00045) = 1.5529$$

PROBLEMS

1. A solution consisting of isobutyl bromide and isobutyl chloride is found to have a refractive index of 1.3931 at 20°C. The refractive indices at 20°C of isobutyl bromide and isobutyl chloride are 1.4368 and 1.3785, respectively. Determine the molar composition (in percent) of the mixture by assuming a linear relation between the refractive index and the molar composition of the mixture.
2. The refractive index of a compound at 16°C is found to be 1.3982. Correct this refractive index to 20°C.

Infrared Spectroscopy

Almost any compound having covalent bonds, whether organic or inorganic, will be found to absorb frequencies of electromagnetic radiation in the infrared region of the spectrum. The infrared region of the electromagnetic spectrum lies at wavelengths longer than those associated with visible light, which includes wavelengths from approximately 400 nm to 800 nm ($1 \text{ nm} = 10^{-9} \text{ m}$), but at wavelengths shorter than those associated with radio waves, which have wavelengths longer than 1 cm. For chemical purposes, we are interested in the *vibrational* portion of the infrared region. This portion includes radiations with wavelengths (λ) between $2.5 \mu\text{m}$ and $15 \mu\text{m}$ ($1 \mu\text{m} = 10^{-6} \text{ m}$). The relation of the infrared region to other regions included in the electromagnetic spectrum is illustrated in Figure 25.1.

As with other types of energy absorption, molecules are excited to a higher energy state when they absorb infrared radiation. The absorption of the infrared radiation is, like other absorption processes, a quantized process. Only selected frequencies (energies) of infrared radiation are absorbed by a molecule. The absorption of infrared radiation corresponds to energy changes on the order of 8–40 kJ/mole (2–10 kcal/mole). Radiation in this energy range corresponds to the range encompassing the stretching and bending vibrational frequencies of the bonds in most covalent molecules. In the absorption process, those frequencies of infrared radiation that match the natural vibrational frequencies of the molecule in question are absorbed, and the energy absorbed increases the *amplitude* of the vibrational motions of the bonds in the molecule.

Most chemists refer to the radiation in the vibrational infrared region of the electromagnetic spectrum by units called **wavenumbers** ($\bar{\nu}$). Wavenumbers are expressed in reciprocal centimeters (cm^{-1}) and are easily computed by taking the reciprocal of the wavelength (λ) expressed in centimeters. This unit has the advantage, for those performing calculations, of being directly proportional to energy. Thus, the vibrational infrared region of the spectrum extends from about 4000 cm^{-1} to 650 cm^{-1} (or wavenumbers).

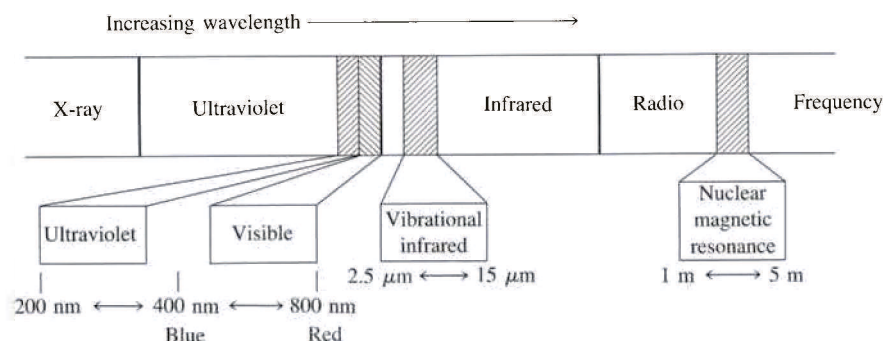


Figure 25.1

A portion of the electromagnetic spectrum showing the relation of vibrational infrared radiation to other types of radiation.

Wavelengths (μm) and wavenumbers (cm^{-1}) can be interconverted by the following relationships:

$$\text{cm}^{-1} = \frac{1}{(\mu\text{m})} \times 10,000$$

$$\mu\text{m} = \frac{1}{(\text{cm})^{-1}} \times 10,000$$

PART A. SAMPLE PREPARATION AND RECORDING THE SPECTRUM

25.1 Introduction

To determine the infrared spectrum of the compound, one must place the compound in a sample holder or cell. In infrared spectroscopy, this immediately poses a problem. Glass, quartz, and plastics absorb strongly throughout the infrared region of the spectrum (any compound with covalent bonds usually absorbs) and cannot be used to construct sample cells. Ionic substances must be used in cell construction. Metal halides (sodium chloride, potassium bromide, silver chloride) are commonly used for this purpose.

Sodium Chloride Cells. Single crystals of sodium chloride are cut and polished to give plates that are transparent throughout the infrared region. These plates are then used to fabricate cells that can be used to hold *liquid* samples. Because sodium chloride is water soluble, samples must be *dry* before a spectrum can be obtained. In general, sodium chloride plates are preferred for most applications involving liquid samples. Potassium bromide plates may also be used in place of sodium chloride.

Silver Chloride Cells. Cells may be constructed of silver chloride. These plates may be used for *liquid* samples that contain small amounts of water, because silver chloride is water-insoluble. However, because water absorbs in the infrared region, as much water as possible should be removed, even when using silver chloride. Silver chloride plates must be stored in the dark. They darken when exposed to light, and they cannot be used with compounds that have an amino functional group. Amines react with silver chloride.

Solid Samples. The easiest way to hold a *solid* sample in place is to dissolve the sample in a volatile organic solvent, place several drops of this solution on a salt plate, and allow the solvent to evaporate. This dry film method can be used only with modern FT-IR spectrometers. The other methods described here can be used with both FT-IR and dispersion spectrometers. A solid sample can also be held in place by making a potassium bromide pellet that contains a small amount of dispersed compound. A solid sample may also be suspended in mineral oil, which absorbs only in specific regions of the infrared spectrum. Another method is to dissolve the solid compound in an appropriate solvent and place the solution between two sodium chloride or silver chloride plates.

ATR Accessory. Modern FT-IR instruments now offer an attenuated total reflectance (ATR) accessory along with the typical transmittance module. The ATR method provides a powerful sampling technique that virtually eliminates sample preparation with both liquids and solids, thus leading to a dramatic improvement in throughput in the teaching laboratory. Although manufacturers offer multiple crystal options, the diamond ATR offers the best option for maximum durability in

the organic teaching laboratory. With the ATR accessory, one simply places a small amount of a liquid or solid directly on the diamond without any previous preparation. The spectrum obtained with an ATR FT-IR is nearly identical to that obtained with an FT-IR operating in the transmittance mode. One may observe some differences in the relative intensities of the peaks, but the peak position in wavenumbers is identical in both modes.¹ ATR FT-IR does not require a clear sample that allows light to pass through the sample such as is common with transmittance instruments. There are some limitations with a diamond ATR instrument. Some materials such as coating on metal and very dark samples do not analyze satisfactorily, but there are few other limitations.

25.2 Liquid Samples—NaCl Plates

The simplest method of preparing the sample, if it is a liquid, is to place a thin layer of the liquid between two sodium chloride plates that have been ground flat and polished. This is the method of choice when you need to determine the infrared spectrum of a pure liquid. A spectrum determined by this method is referred to as a **neat** spectrum. No solvent is used. The polished plates are expensive because they are cut from a large, single crystal of sodium chloride. Salt plates break easily, and they are water soluble.

Preparing the Sample. Obtain two sodium chloride plates and a holder from the desiccator where they are stored. Moisture from fingers will mar and occlude the polished surfaces. Samples that contain water will destroy the plates.

NOTE: The plates should be touched only on their edges. Be certain to use a sample that is dry or free from water.

Add 1 or 2 drops of the liquid to the surface of one plate and then place the second plate on top.² The pressure of this second plate causes the liquid to spread out and form a thin capillary film between the two plates. As shown in Figure 25.2, set the plates between the bolts in a holder and place the metal ring carefully on the salt plates. Use the hex nuts to hold the salt plates in place.

NOTE: Do not overtighten the nuts or the salt plates will cleave or split.

Tighten the nuts firmly, but do not use any force to turn them. Spin them with the fingers until they stop; then turn them just another fraction of a full turn, and they will be tight enough. If the nuts have been tightened carefully, you should observe a *transparent film of sample* (a uniform wetting of the surface). If a thin film has not been obtained, either loosen one or more of the hex nuts and adjust them so that a uniform film is obtained or add more sample.

The thickness of the film obtained between the two plates is a function of two factors: (1) the amount of liquid placed on the first plate (1 drop, 2 drops, and so on), and (2) the pressure used to hold the plates together. If more than 1 or 2 drops of liquid have been used, the amount will probably be too much, and the resulting

¹Shuttlefield, J. D.; Grassian, V. H. ATR-FTIR in the Undergraduate Chemistry Laboratory, Part 1: "Fundamentals and Examples". *J. Chem. Educ.* 85, (2008): 279–281.

² Use a Pasteur pipette or a short length of microcapillary tubing. If you use the microcapillary tubing, it can be filled by touching it into the liquid sample. When you touch it (lightly) to the salt plate, it will empty. Be careful not to scratch the plate.

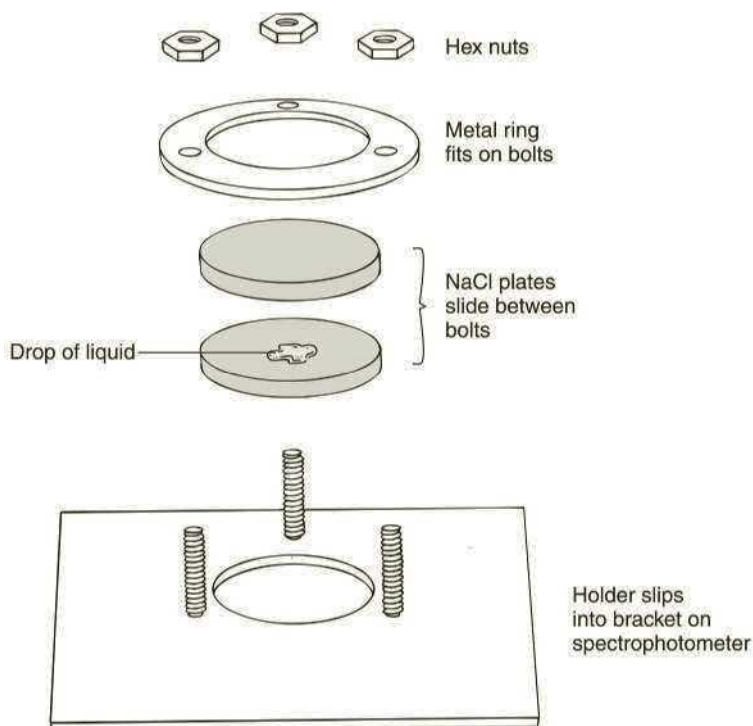


Figure 25.2
Salt plates and holder.

spectrum will show strong absorptions that are off the scale of the chart paper. Only enough liquid to wet both surfaces is needed.

If the sample has a very low viscosity, the capillary film may be too thin to produce a good spectrum. Another problem you may find is that the liquid is so volatile that the sample evaporates before the spectrum can be determined. In these cases, you may need to use the silver chloride plates discussed in Section 25.3 or a solution cell described in Section 25.6. Often, you can obtain a reasonable spectrum by assembling the cell quickly and running the spectrum before the sample runs out of the salt plates or evaporates.

Determining the Infrared Spectrum. Slide the holder into the slot in the sample beam of the spectrophotometer. Determine the spectrum according to the instructions provided by your instructor. In some cases, your instructor may ask you to calibrate your spectrum. If this is the case, refer to Section 25.8.

Cleaning and Storing the Salt Plates. Once the spectrum has been determined, demount the holder and rinse the salt plates with methylene chloride (or *dry* acetone). (Keep the plates away from water!) Use a soft tissue, moistened with the solvent, to wipe the plates. If some of your compound remains on the plates, you may observe a shiny surface. Continue to clean the plates with solvent until no more compound remains on the surfaces of the plates.

CAUTION

Avoid direct contact with methylene chloride. Return the salt plates and holder to the desiccator for storage.



25.3 Liquid Samples—AgCl Plates

The minicell shown in Figure 25.3 may also be used with liquids.³ The cell assembly consists of a two-piece threaded body, an O-ring, and two silver chloride plates. The plates are flat on one side, and there is a circular depression (0.025 mm or 0.10 mm deep) on the other side of the plate. An advantage of using silver chloride plates is that they may be used with wet samples or solutions. A disadvantage is that silver chloride darkens when exposed to light for extended periods. Silver chloride plates also scratch more easily than salt plates and react with amines.

Preparing the Sample. Silver chloride plates should be handled in the same way as salt plates. Unfortunately, they are smaller and thinner (about like a contact lens) than salt plates, and care must be taken not to lose them! Remove them from the light-tight container with care. It is difficult to tell which side of the plate has the slight circular depression. Your instructor may have etched a letter on each plate to indicate which side is the flat one. To determine the infrared spectrum of a pure liquid (neat spectrum), select the flat side of each silver chloride plate. Insert the O-ring into the cell body as shown in Figure 25.3, place the plate into the cell body with the flat surface up, and add 1 drop or less of liquid to the plate.

NOTE: Do not use amines with AgCl plates.

Place the second plate on top of the first with the flat side down. The orientation of the silver chloride plates is shown in Figure 25.4A. This arrangement is used to obtain a capillary film of your sample. Screw the top of the minicell into the body of

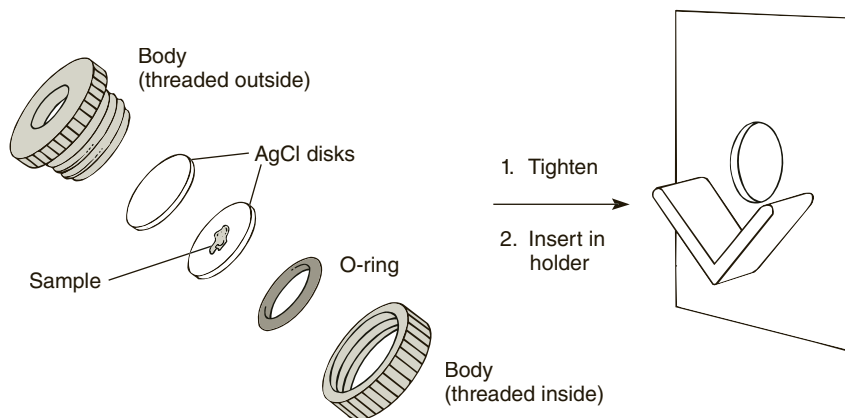


Figure 25.3
AgCl minicell liquid sample cell and V-mount holder.

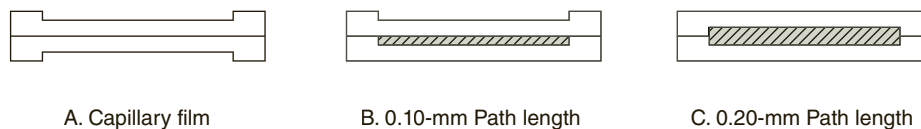


Figure 25.4
Path-length variations for AgCl plates.

³ The Wilks Mini-Cell liquid sample holder is available from the Foxboro Company, 151 Woodward Avenue, South Norwalk, CT 06856. We recommend the AgCl cell windows with 0.10-mm depression rather than the 0.025-mm depression.

the cell so that the silver chloride plates are held firmly together. A tight seal forms because AgCl deforms under pressure.

Other combinations may be used with these plates. For example, you may vary the sample path length by using the orientations shown in Figures 25.4B and C. If you add your sample and the 0.10-mm depression of one plate and cover it with the flat side of the other one, you obtain a path length of 0.10 mm (see Figure 25.4B). This arrangement is useful for analyzing volatile or low-viscosity liquids. Placement of the two plates with their depressions toward each other gives a path length of 0.20 mm (see Figure 25.4C). This orientation may be used for a solution of a solid (or liquid) in carbon tetrachloride (see Section 25.6B).

Determining the Spectrum. Slide the V-mount holder shown in Figure 25.3 into the slot on the infrared spectrophotometer. Set the cell assembly in the V-mount holder and determine the infrared spectrum of the liquid.

Cleaning and Storing the AgCl Plates. Once the spectrum has been determined, the cell assembly holder should be demounted and the AgCl plates rinsed with methylene chloride or acetone. Do not use tissue to wipe the plates, because they scratch easily. AgCl plates are light sensitive. Store the plates in a light-tight container.

25.4 Solid Samples—Dry Film

A simple method for determining the infrared spectrum of a solid sample is the dry film method. This method is easier than the other methods described here, it does not require any specialized equipment, and the spectra are excellent.⁴ The disadvantage is that the dry film method can be used only with modern FT-IR spectrometers.

To use this method, place about 5 mg of your solid sample in a small, clean test tube. Add about 5 drops of methylene chloride (or diethyl ether, pentane, or dry acetone), and stir the mixture to dissolve the solid. Using a Pasteur pipette (not a capillary tube), place several drops of the solution on the face of a salt plate. Allow the solvent to evaporate; a uniform deposit of your product will remain as a dry film coating the salt plate. Mount the salt plate on a V-shaped holder in the infrared beam. Note that only one salt plate is used; the second salt plate is not used to cover the first. Once the salt plate is positioned properly, you may determine the spectrum in the normal manner. With this method, it is *very important* that you clean your material off the salt plate. When you are finished, use methylene chloride or dry acetone to clean the salt plate.

25.5 Solid Samples—Other Methods

The methods described in this section can be used with both FT-IR and dispersion spectrometers.

A. Attenuated Total Reflectance (ATR)

Preparation of solid samples for analysis by infrared analysis has tended to be much more labor intensive than analysis of liquids. The traditional methods involve preparation of a KBr pellet or a Nujol mull described in sections 25.5 B and 25.5 C. The use of an attenuated total reflectance accessory with an FT-IR instrument has dramatically improved the preparation time previously required. One simply places a small amount of a solid on the instrument. For these laboratories equipped with this accessory, it is strongly recommended for analysis of both solids and liquids (Section 25.1).

⁴ Feist, P. L. Sampling Techniques for Organic Solids in IR Spectroscopy: Thin Solid Films as the Method of Choice in Teaching Laboratories. *Journal of Chemical Education*, 78 (2001): 351.

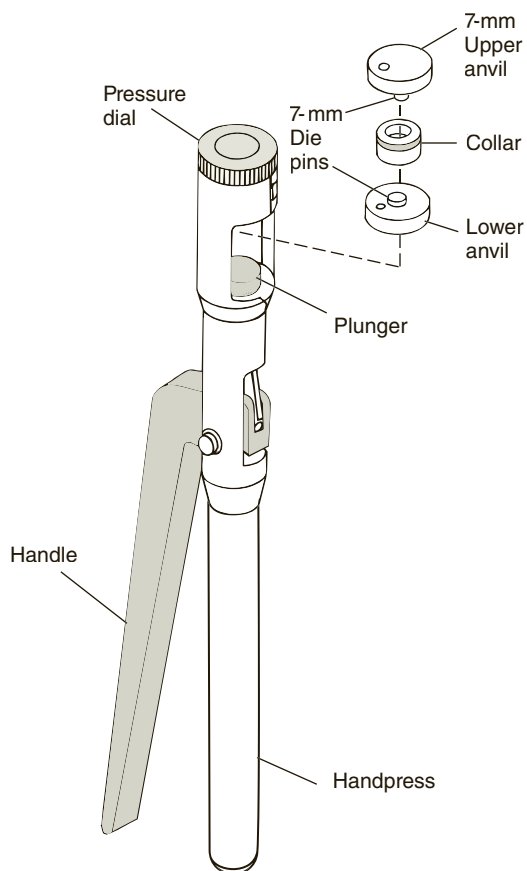


Figure 25.5
Making a KBr pellet with a handpress.

B. KBr Pellets

One method of preparing a solid sample is to make a **potassium bromide (KBr) pellet**. When KBr is placed under pressure, it melts, flows, and seals the sample into a solid solution, or matrix. Because potassium bromide does not absorb in the infrared spectrum, a spectrum can be obtained on a sample without interference.

Preparing the Sample. Remove the agate mortar and pestle from the desiccator for use in preparing the sample. (Take care of them; they are expensive.) Grind 1 mg (0.001 g) of the solid sample for 1 minute in the agate mortar. At this point, the particle size will become so small that the surface of the solid appears shiny. Add 80 mg (0.080 g) of *powdered* KBr and grind the mixture for about 30 seconds with the pestle. Scrape the mixture into the middle with a spatula and grind the mixture again for about 15 seconds. This grinding operation helps to mix the sample thoroughly with the KBr. You should work as rapidly as possible because KBr absorbs water. The sample and KBr must be finely ground, or the mixture will scatter the infrared radiation excessively. Using your spatula, heap the mixture in the center of the mortar. Return the bottle of potassium bromide to the desiccator where it is stored when it is not in use.

The sample and potassium bromide should be weighed on an analytical balance the first few times that a pellet is prepared. After some experience, you can estimate these quantities quite accurately by eye.

Making a Pellet Using a KBr Handpress. Two methods are commonly used to prepare KBr pellets. The first method uses the handpress apparatus shown in Figure 25.5.⁵ Remove the die set from the storage container. Take extreme care to avoid scratching the polished surfaces of the die set. Place the anvil with the shorter die pin (lower anvil in Figure 25.5) on a bench. Slip the collar over the pin. Remove about one fourth of your KBr mixture with a spatula and transfer it into the collar. The powder may not cover the head of the pin completely, but do not be concerned about this. Place the anvil with the longer die pin into the collar so that the die pin comes into contact with the sample. Never press the die set unless it contains a sample.

Lift the die set carefully by holding onto the lower anvil so that the collar stays in place. If you are careless with this operation, the collar may move enough to allow the powder to escape. Open the handle of the handpress slightly, tilt the press back a bit, and insert the die set into the press. Make sure that the die set is seated against the side wall of the chamber. Close the handle. It is imperative that the die set be seated against the side wall of the chamber so that the die is centered in the chamber. Pressing the die in an off-centered position can bend the anvil pins.

With the handle in the closed position, rotate the pressure dial so that the upper ram of the handpress just touches the upper anvil of the die assembly. Tilt the unit back so that the die set does not fall out of the handpress. Open the handle and rotate the pressure dial clockwise about one-half turn. Slowly compress the KBr

⁵The KBr Quick Press unit is available from Wilmad Glass Company, Inc., Route 40 and Oak Road, Buena, NJ 08310.

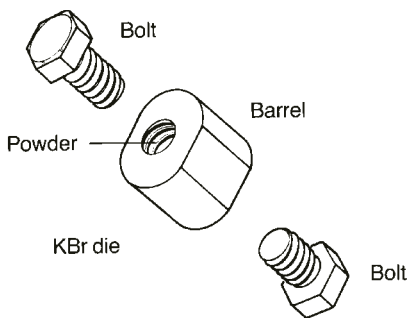


Figure 25.6
Making a KBr pellet with a minipress.

mixture by closing the handle. The pressure should be no greater than that exerted by a very firm handshake. Do not apply excessive pressure or the dies may be damaged. If in doubt, rotate the pressure dial counterclockwise to lower the pressure. If the handle closes too easily, open the handle, rotate the pressure dial clockwise, and compress the sample again. Compress the sample for about 60 seconds.

After this time, tilt the unit back so that the die set does not fall out of the handpress. Open the handle and carefully remove the die set from the unit. Turn the pressure dial counterclockwise about one full turn. Pull the die set apart and inspect the KBr pellet. Ideally, the pellet should appear clear like a piece of glass, but usually it will be translucent or somewhat opaque. There may be some cracks or holes in the pellet. The pellet will produce a good spectrum, even with imperfections, as long as light can travel through the pellet. Clean the dies using the procedure outlined below, in “Cleaning and Storing the Equipment.”

Making a Pellet with a KBr Minipress. The second method of preparing a pellet uses the minipress apparatus shown in Figure 25.6. Obtain a ground KBr mixture as described in “Preparing the Sample” and transfer a portion of the finely ground powder (usually not more than half) into a die that compresses it into a translucent pellet. As shown in Figure 25.6, the die consists of two stainless steel bolts and a threaded barrel. The bolts have their ends ground flat. To use this die, screw one of the bolts into the barrel, but not all the way; leave one or two turns. Carefully add the powder with a spatula into the open end of the partly assembled die and tap it lightly on the benchtop to give an even layer on the face of the bolt. While keeping the barrel upright, carefully screw the second bolt into the barrel until it is finger tight. Insert the head of the bottom bolt into the hexagonal hole in a plate bolted to the benchtop. This plate keeps the head of one bolt from turning. The top bolt is tightened with a torque wrench to compress the KBr mixture. Continue to turn the torque wrench until you hear a loud click (the ratchet mechanism makes softer clicks) or until you reach the appropriate torque value (20 ft-lb). If you tighten the bolt beyond this point, you may twist the head off one of the bolts. Leave the die under pressure for about 60 seconds; then reverse the ratchet on the torque wrench or pull the torque wrench in the opposite direction to open the assembly. When the two bolts are loose, hold the barrel horizontally and carefully remove the two bolts. You should observe a clear or translucent KBr pellet in the center of the barrel. Even if the pellet is not totally transparent, you should be able to obtain a satisfactory spectrum as long as light passes through the pellet.

Determining the Infrared Spectrum. To obtain the spectrum, slide the holder appropriate for the type of die that you are using into the slot on the infrared spectrophotometer. Set the die containing the pellet in the holder so that the sample is centered in the optical path. Obtain the infrared spectrum. If you are using a double-beam instrument, you may be able to compensate (at least partially) for a marginal pellet by placing a wire screen or attenuator in the reference beam, thereby balancing the lowered transmittance of the pellet. An FT-IR instrument will automatically deal with the low intensity if you select the “autoscale” option.

Problems with an Unsatisfactory Pellet. If the pellet is unsatisfactory (too cloudy to pass light), one of several things may have been wrong:

1. The KBr mixture may not have been ground finely enough, and the particle size may be too big. The large particle size creates too much light scattering.
2. The sample may not be dry.

3. Too much sample may have been used for the amount of KBr taken.
4. The pellet may be too thick; that is, too much of the powdered mixture was put into the die.
5. The KBr may have been “wet” or have acquired moisture from the air while the mixture was being ground in the mortar.
6. The sample may have a low melting point. Low-melting solids not only are difficult to dry, but also melt under pressure. You may need to dissolve the compound in a solvent and run the spectrum in solution (see Section 25.6).

Cleaning and Storing the Equipment. After you have determined the spectrum, punch the pellet out of the die with a wooden applicator stick (a spatula should not be used as it may scratch the dies). Remember that the polished faces of the die set must not be scratched, or they become useless. After the pellet has been punched out, wash all parts of the die set or minipress with warm water. Then rinse the parts with acetone and dry them using a Kimwipe. Check with your instructor to see if there are additional instructions for cleaning the die set. Return the dies to the storage container. Wash the mortar and pestle with water, dry them carefully with paper towels, and return them to the desiccator. Return the KBr powder to its desiccator.

C. Nujol Mulls

If an adequate KBr pellet cannot be obtained or if the solid is insoluble in a suitable solvent, the spectrum of a solid may be determined as a **Nujol mull**. In this method, finely grind about 5 mg of the solid sample in an agate mortar with a pestle. Then add 1 or 2 drops of Nujol mineral oil (white) and grind the mixture to a very fine dispersion. The solid is not dissolved in the Nujol; it is actually a suspension. This mull is then placed between two salt plates using a rubber policeman. Mount the salt plates in the holder in the same way as for liquid samples (see Section 25.2).

Nujol is a mixture of high-molecular-weight hydrocarbons. Hence, it has absorptions in the C—H stretch and CH₂ and CH₃ bending regions of the spectrum (Figure 25.7). Clearly, if Nujol is used, no information can be obtained in these portions of the spectrum. In interpreting the spectrum, you must ignore these Nujol peaks. It is important to label the spectrum immediately after it was determined,

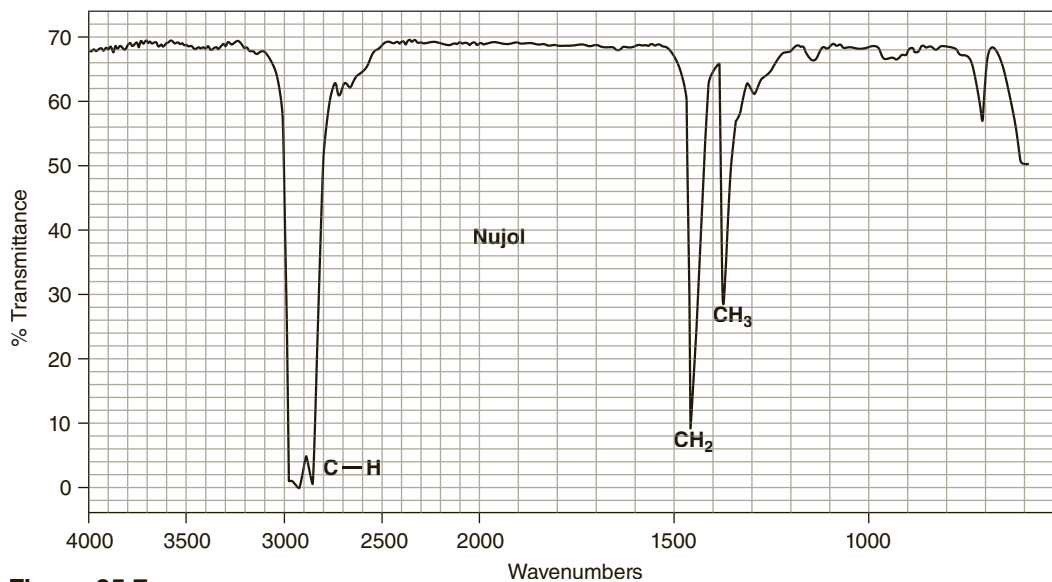


Figure 25.7
Infrared spectrum of Nujol (mineral oil).

noting that it was determined as a Nujol mull. Otherwise, you might forget that the C—H peaks belong to Nujol and not to the dispersed solid.

25.6 Solid Samples—Solution Spectra

A. Method A—Solution Between Salt (NaCl) Plates

For substances that are soluble in carbon tetrachloride, a quick and easy method for determining the spectra of solids is available. Dissolve as much solid as possible in 0.1 mL of carbon tetrachloride. Place 1 or 2 drops of the solution between sodium chloride plates in precisely the same manner as used for pure liquids (see Section 25.2). The spectrum is determined as described for pure liquids using salt plates (see Section 25.2). You should work as quickly as possible. If there is a delay, the solvent will evaporate from between the plates before the spectrum is recorded. Because the spectrum contains the absorptions of the solute superimposed on the absorptions of carbon tetrachloride, it is important to remember that any absorption that appears near 800 cm^{-1} may be due to the stretching of the C—Cl bond of the solvent. Information contained to the right of about 900 cm^{-1} is not usable in this method. There are no other interfering bands for this solvent (see Figure 25.8), and any other absorptions can be attributed to your sample. Chloroform solutions should not be studied by this method because the solvent has too many interfering absorptions (see Figure 25.9).

CAUTION



Carbon tetrachloride is a hazardous solvent. Work under the hood!

Carbon tetrachloride, besides being toxic, is suspected of being a carcinogen. Despite the health problems associated with its use, there is no suitable alternative solvent for infrared spectroscopy. Other solvents have too many interfering infrared absorption bands. Handle carbon tetrachloride carefully to minimize the adverse health effects. The spectroscopic-grade carbon tetrachloride should be stored in a glass-stoppered bottle in a hood. A Pasteur pipette should be attached to the bottle, possibly by storing it in a test tube taped to the side of the bottle. All sample preparation should be conducted in a hood. Rubber or plastic gloves should be worn.

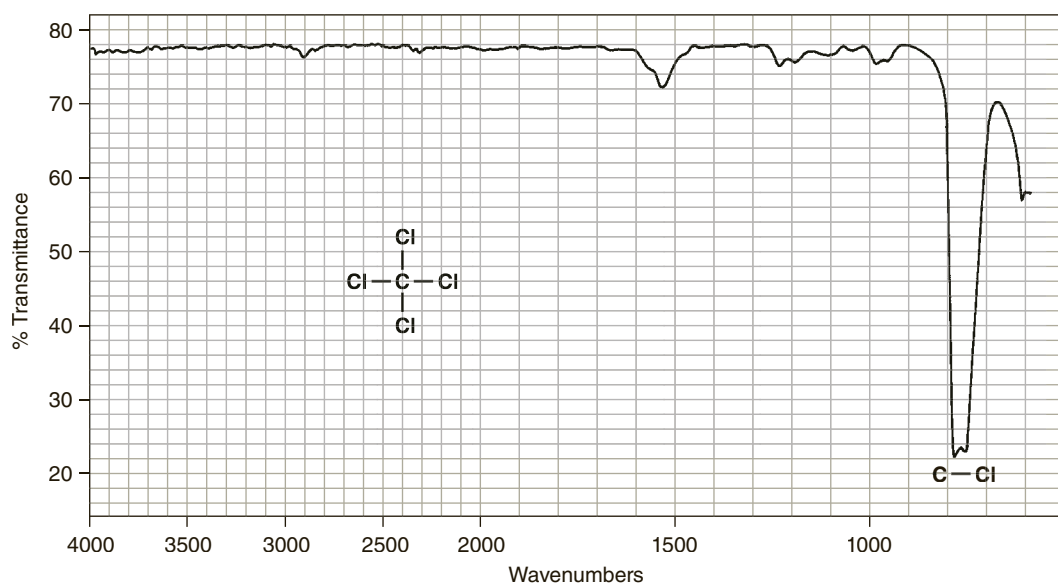


Figure 25.8
Infrared spectrum of carbon tetrachloride.

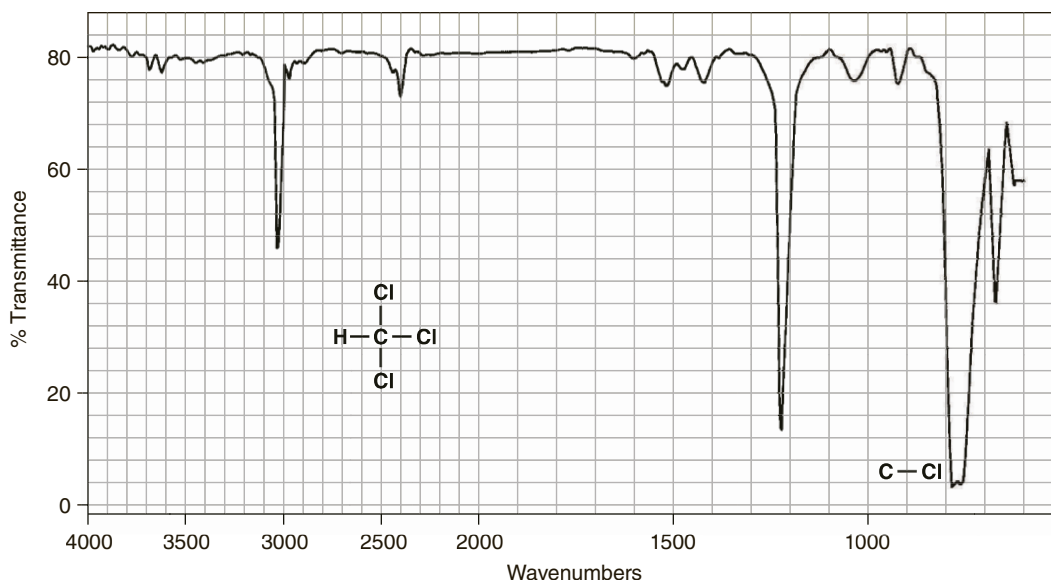


Figure 25.9
Infrared spectrum of chloroform.

The cells should also be cleaned in the hood. All carbon tetrachloride used in preparing samples should be disposed of in an appropriately marked waste container.

B. Method B—AgCl Minicell

The AgCl minicell described in Section 25.3 may be used to determine the infrared spectrum of a solid dissolved in carbon tetrachloride. Prepare a 5–10% solution (5–10 mg in 0.1 mL) in carbon tetrachloride. If it is not possible to prepare a solution of this concentration because of low solubility, dissolve as much solid as possible in the solvent. Following the instructions given in Section 25.3, position the AgCl plates as shown in Figure 25.4C to obtain the maximum possible path length of 0.20 mm. When the cell is tightened firmly, the cell will not leak.

As indicated in method A, the spectrum will contain the absorptions of the dissolved solid superimposed on the absorptions of carbon tetrachloride. A strong absorption appears near 800 cm^{-1} for the C—Cl stretch in the solvent. No useful information may be obtained for the sample to the right of about 900 cm^{-1} , but other bands that appear in the spectrum will belong to your sample. Read the safety material provided in method A. Carbon tetrachloride is toxic, and it should be used under a hood.

NOTE: Care should be taken in cleaning the AgCl plates. Because AgCl plates scratch easily, they should not be wiped with tissue. Rinse them with methylene chloride and keep them in a dark place. Amines will destroy the plates.

C. Method C—Solution Cells (NaCl)

The spectra of solids may also be determined in a type of permanent sample cell called a **solution cell**. (The infrared spectra of liquids may also be determined in this cell.) The solution cell, shown in Figure 25.10, is made from two salt plates, mounted with a Teflon spacer between them to control the thickness of the sample. The top sodium chloride plate has two holes drilled in it so that the sample can be introduced into the cavity between the two plates. These holes are extended

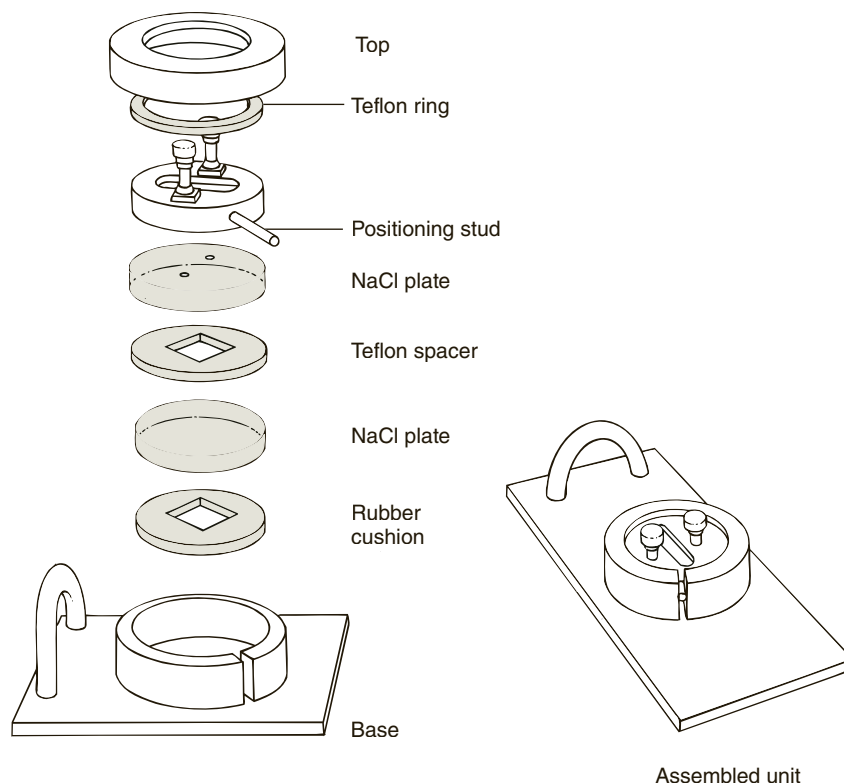


Figure 25.10
A solution cell.

through the face plate by two tubular extensions designed to hold Teflon plugs, which seal the internal chamber and prevent evaporation. The tubular extensions are tapered so that a syringe body (Luer lock without a needle) will fit snugly into them from the outside. The cells are thus filled from a syringe; usually, they are held upright and filled from the bottom entrance port.

These cells are expensive, and you should try either method A or B before using solution cells. If you do need them, obtain your instructor's permission and receive instruction before using the cells. The cells are purchased in matched pairs, with identical path lengths. Dissolve a solid in a suitable solvent, usually carbon tetrachloride, and add the solution to one of the cells (**sample cell**) as described in the previous paragraph. The pure solvent, identical to that used to dissolve the solid, is placed in the other cell (**reference cell**). The spectrum of the solvent is subtracted from the spectrum of the solution (not always completely), and a spectrum of the solute is thus provided. For the solvent compensation to be as exact as possible and to avoid contamination of the reference cell, it is essential that one cell be used as a reference and that the other cell be used as a sample cell without ever being interchanged. After the spectrum is determined, it is important to clean the cells by flushing them with clean solvent. They should be dried by passing dry air through the cell.

Solvents most often used in determining infrared spectra are carbon tetrachloride (see Figure 25.8), chloroform (see Figure 25.9), and carbon disulfide (see Figure 25.11). A 5–10% solution of solid in one of these solvents usually gives a good spectrum. Carbon tetrachloride and chloroform are suspected carcinogens; however, because there are no suitable alternative solvents, these compounds must be used in infrared spectroscopy. The procedure outlined above for carbon tetrachloride should be followed. This procedure serves equally well for chloroform.

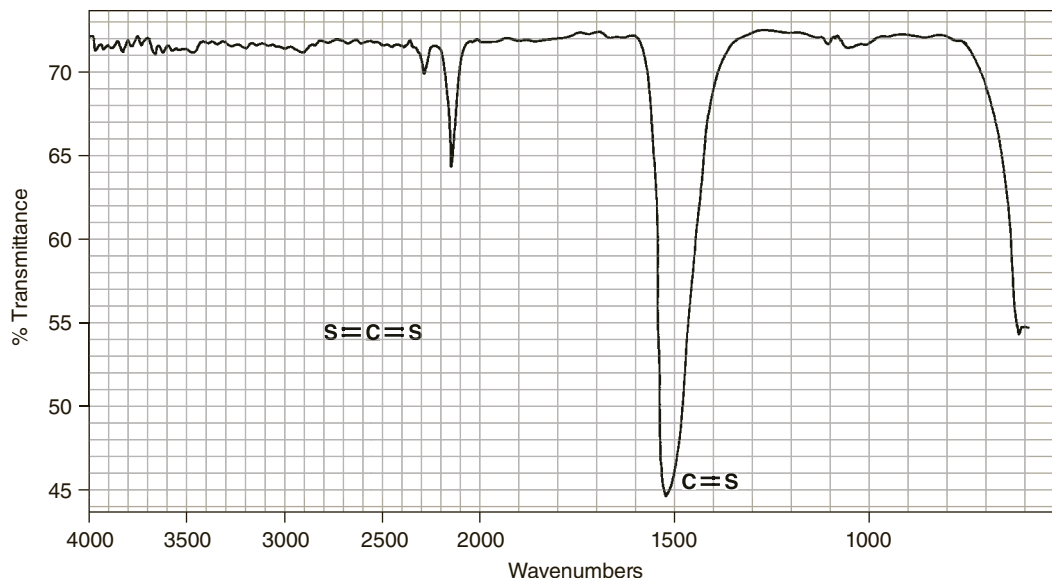


Figure 25.11
Infrared spectrum of carbon disulfide.

NOTE: Before you use the solution cells, you must obtain the instructor's permission and instruction on how to fill and clean the cells.

25.7 Recording the Spectrum

The instructor will describe how to operate the infrared spectrophotometer, because the controls vary considerably, depending on the manufacturer, model of the instrument, and type. For example, some instruments involve pushing only a few buttons, whereas others use a more complicated computer interface system.

In all cases, it is important that the sample, the solvent, the type of cell or method used, and any other pertinent information be written on the spectrum immediately after the determination. This information may be important, and it is easily forgotten if not recorded. You may also need to calibrate the instrument (see Section 25.8).

25.8 Calibration

For some instruments, the frequency scale of the spectrum must be calibrated so that you know the position of each absorption peak precisely. You can recalibrate by recording a very small portion of the spectrum of polystyrene over the spectrum of your sample. The complete spectrum of polystyrene is shown in Figure 25.12. The most important of these peaks is at 1603 cm^{-1} ; other useful peaks are at 2850 cm^{-1} and 906 cm^{-1} . After you record the spectrum of your sample, substitute a thin film of polystyrene for the sample cell and record the tips (not the entire spectrum) of the most important peaks over the sample spectrum.

It is always a good idea to calibrate a spectrum when the instrument uses chart paper with a preprinted scale. It is difficult to align the paper properly so that the scale matches the absorption lines precisely. You often need to know the precise values for certain functional groups (for example, the carbonyl group). Calibration is essential in these cases.

With computer-interfaced instruments, the instrument does not need to be calibrated. With this type of instrument, the spectrum and scale are printed on blank

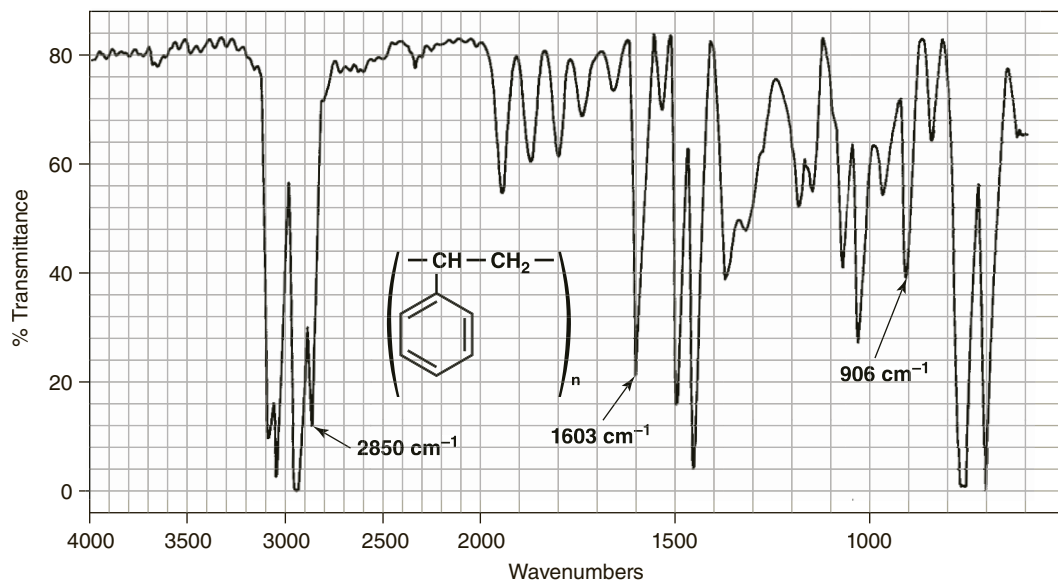


Figure 25.12
Infrared spectrum of polystyrene (thin film).


paper at the same time. The instrument has an internal calibration that ensures that the positions of the absorptions are known precisely and that they are placed at the proper positions on the scale. With this type of instrument, it is often possible to print a list of the locations of the major peaks as well as to obtain the complete spectrum of your compound.

PART B. INFRARED SPECTROSCOPY

25.9 Uses of the Infrared Spectrum

Because every type of bond has a different natural frequency of vibration and because the same type of bond in two different compounds is in a slightly different environment, no two molecules of different structure have exactly the same infrared absorption pattern, or **infrared spectrum**. Although some of the frequencies absorbed in the two cases might be the same, in no case of two different molecules will their infrared spectra (the patterns of absorption) be identical. Thus, the infrared spectrum can be used to identify molecules much as a fingerprint can be used to identify people. Comparing the infrared spectra of two substances thought to be identical will establish whether or not they are in fact identical. If the infrared spectra of two substances coincide peak for peak (absorption for absorption), in most cases, the substances are identical.


A second and more important use of the infrared spectrum is that it gives structural information about a molecule. The absorptions of each type of bond (N—H, C—H, O—H, C—X, C=O, C—O, C—C, C=C, C≡C, C≡N, and so on) are regularly found only in certain small portions of the vibrational infrared region. A small range of absorption can be defined for each type of bond. Outside this range, absorptions will normally be due to some other type of bond. Thus, for instance, any absorption in the range $3000 \pm 150 \text{ cm}^{-1}$ will almost always be due to the presence of a CH bond in the molecule; an absorption in the range $1700 \pm 100 \text{ cm}^{-1}$ will normally be due to the presence of a C=O bond (carbonyl group) in the molecule. The same type of range applies to each type of bond. The way these are spread out over the vibrational infrared is illustrated schematically in Figure 25.13. It is a good idea to remember this general scheme for future convenience.

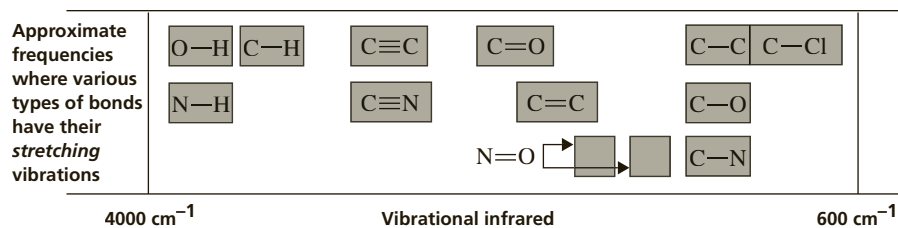**Figure 25.13**

Approximate regions in which various common types of bonds absorb. (Bending, twisting, and other types of bond vibration have been omitted for clarity.)


25.10 Modes of Vibration

The simplest types, or **modes**, of vibrational motion in a molecule that are **infrared active**, that is, give rise to absorptions, are the stretching and bending modes.

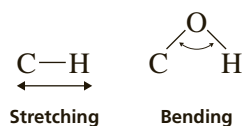

Other, more complex types of stretching and bending are also active, however. To introduce several words of terminology, the normal modes of vibration for a methylene group are shown below.

In any group of three or more atoms—at least two of which are identical—there are *two* modes of stretching or bending: the symmetric mode and asymmetric mode. Examples of such groupings are $-\text{CH}_3$, $-\text{CH}_2-$, $-\text{NO}_2$, $-\text{NH}_2$, and anhydrides $(\text{CO})_2\text{O}$. For the anhydride, owing to asymmetric and symmetric modes of stretch, this functional group gives *two* absorptions in the $\text{C}=\text{O}$ region. A similar phenomenon is seen for amino groups, where primary amines usually have *two* absorptions in the NH stretch region, whereas secondary amines R_2NH have only one absorption peak. Amides show similar bands. There are two strong $\text{N}=\text{O}$ stretch peaks for a nitro group, which are caused by asymmetric and symmetric stretching modes.

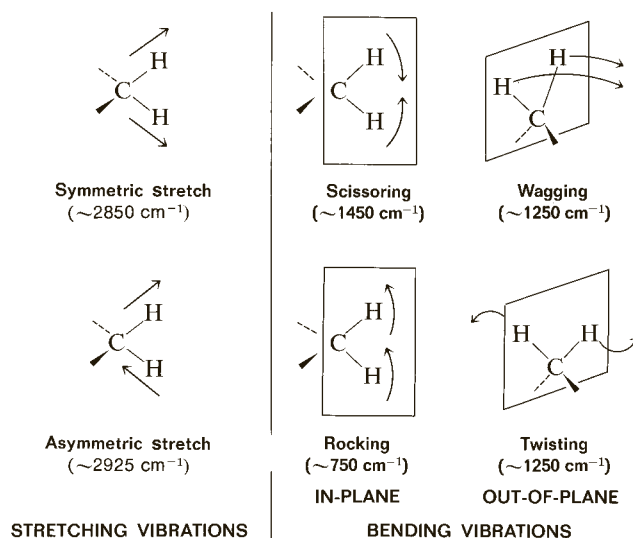

25.11 What to Look for in Examining Infrared Spectra

The instrument that determines the absorption spectrum for a compound is called an **infrared spectrophotometer**. The spectrophotometer determines the relative strengths and positions of all the absorptions in the infrared region and

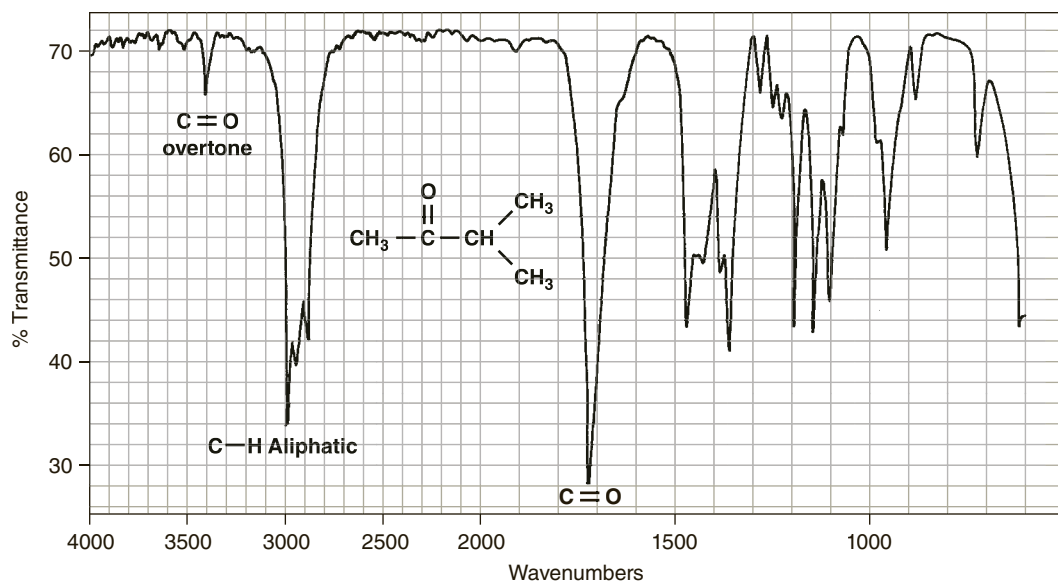


Figure 25.14
Infrared spectrum of methyl isopropyl ketone (neat liquid, salt plates).

plots this information on a piece of paper. This plot of absorption intensity versus wavenumber or wavelength is referred to as the **infrared spectrum** of the compound. A typical infrared spectrum, that of methyl isopropyl ketone, is shown in Figure 25.14.

The strong absorption in the middle of the spectrum corresponds to C=O, the carbonyl group. Note that the C=O peak is quite intense. In addition to the characteristic position of absorption, the **shape** and **intensity** of this peak are also unique to the C=O bond. This is true for almost every type of absorption peak; both shape and intensity characteristics can be described, and these characteristics often make it possible to distinguish the peak in a confusing situation. For instance, to some extent both C=O and C=C bonds absorb in the same region of the infrared spectrum:

$$\text{C=O } 1850\text{--}1630 \text{ cm}^{-1}$$

$$\text{C=C } 1680\text{--}1620 \text{ cm}^{-1}$$

However, the C=O bond is a strong absorber, whereas the C=C bond generally absorbs only weakly. Hence, a trained observer would not normally interpret a strong peak at 1670 cm^{-1} to be a carbon-carbon double bond or a weak absorption at this frequency to be due to a carbonyl group.

In Figure 25.14, notice that the very strong C=O peak at 1710 cm^{-1} leads to a very weak overtone at twice the C=O frequency (3420 cm^{-1}). As is often the case, very strong bands may yield a weaker peak at twice the frequency of the main band.

The shape of a peak often gives a clue to its identity as well. Thus, although the NH and OH regions of the infrared overlap,

$$\text{OH } 3650\text{--}3200 \text{ cm}^{-1}$$

$$\text{NH } 3500\text{--}3300 \text{ cm}^{-1}$$

NH usually gives a **sharp** absorption peak (absorbs a very narrow range of frequencies), and OH, when it is in the NH region, usually gives a **broad** absorption

peak. Primary amines give *two* absorptions in this region, whereas alcohols give only one.

Therefore, while you are studying the sample spectra in the pages that follow, you should also notice shapes and intensities. They are as important as the frequency at which an absorption occurs, and you must train your eye to recognize these features. In the literature of organic chemistry, you will often find absorptions referred to as strong (s), medium (m), weak (w), broad, or sharp. The author is trying to convey some idea of what the peak looks like without actually drawing the spectrum. Although the intensity of an absorption often provides useful information about the identity of a peak, be aware that the relative intensities of all of the peaks in the spectrum are dependent on the amount of sample that is used and the sensitivity setting of the instrument. Therefore, the *actual* intensity of a particular peak may vary from spectrum to spectrum, and you must pay attention to *relative* intensities.

25.12 Correlation Charts and Tables

To extract structural information from infrared spectra, you must know the frequencies or wavelengths at which various functional groups absorb. Infrared **correlation tables** present as much information as is known about where the various functional groups absorb. The books listed at the end of this chapter present extensive lists of correlation tables. Sometimes, the absorption information is given in a chart, called a **correlation chart**. A simplified correlation table is given in Table 25.1.

Although you may think assimilating the mass of data in Table 25.1 will be difficult, it is not if you make a modest start and then gradually increase your familiarity with the data. An ability to interpret the fine details of an infrared spectrum will follow. This is most easily accomplished by first establishing the broad visual patterns of Figure 25.13 firmly in mind. Then, as a second step, a “typical absorption value” can be memorized for each of the functional groups in this pattern. This value will be a single number that can be used as a pivot value for the memory. For instance, start with a simple aliphatic ketone as a model for all typical carbonyl compounds. The typical aliphatic ketone has a carbonyl absorption of $1715 \pm 10 \text{ cm}^{-1}$. Without worrying about the variation, memorize 1715 cm^{-1} as the base value for carbonyl absorption. Then learn the extent of the carbonyl range and the visual pattern of how the different kinds of carbonyl groups are arranged throughout this region. See, for instance, Figure 25.27, which gives typical values for carbonyl compounds. Also learn how factors such as ring size (when the functional group is contained in a ring) and conjugation affect the base values (that is, in which direction the values are shifted). Learn the trends—always remembering the base value (1715 cm^{-1}). It might prove useful as a beginning to memorize the base values in Table 25.2 for this approach. Notice that there are only eight values.

25.13 Analyzing a Spectrum (or What You Can Tell at a Glance)

In analyzing the spectrum of an unknown, concentrate first on establishing the presence (or absence) of a few major functional groups. The most conspicuous peaks are C=O, O—H, N—H, C—O, C=C, C≡C, C≡N, and NO₂. If they are present, they give immediate structural information. Do not try to analyze in detail the CH absorptions near 3000 cm^{-1} ; almost all compounds have these absorptions. Do not worry about subtleties of the exact type of environment in which the functional group is found. A checklist of the important gross features follows:

1. Is a carbonyl group present?

Table 25.1 A Simplified Correlation Table

	<i>Type of Vibration</i>		<i>Frequency (cm⁻¹)</i>	<i>Intensity^a</i>
C—H	Alkanes	(stretch)	3000–2850	s
	—CH ₃	(bend)	1450 and 1375	m
	—CH ₂ ⁻	(bend)	1465	m
	Alkenes	(stretch)	3100–3000	m
		(bend)	1700–1000	s
	Aromatics	(stretch)	3150–3050	s
		(out-of-plane bend)	1000–700	s
	Alkyne	(stretch)	ca. 3300	s
Aldehyde			2900–2800	w
			2800–2700	w
C—C	Alkane	Not interpretatively useful		
C=C	Alkene		1680–1600	m–w
	Aromatic		1600–1400	m–w
C≡C	Alkyne		2250–2100	m–w
C=O	Aldehyde		1740–1720	s
	Ketone (acyclic)		1725–1705	s
	Carboxylic acid		1725–1700	s
	Ester		1750–1730	s
	Amide		1700–1640	s
	Anhydride		ca. 1810	s
			ca. 1760	s
C—O	Alcohols, ethers, esters, carboxylic acids		1300–1000	s
O—H	Alcohol, phenols			
	Free		3650–3600	m
	H-Bonded		3400–3200	m
	Carboxylic acids		3300–2500	m
N—H	Primary and secondary amines		ca. 3500	m
C≡N	Nitriles		2260–2240	m
N=O	Nitro (R—NO ₂)		1600–1500	s
			1400–1300	s
C—X	Fluoride		1400–1000	s
	Chloride		800–600	s
	Bromide, iodide		<600	s

^as, strong; m, medium; w, weak.

Table 25.2 Base Values for Absorptions of Bonds

O—H	3400 cm ⁻¹	C≡C	2150 cm ⁻¹
N—H	3500 cm ⁻¹	C=O	1715 cm ⁻¹
C—H	3000 cm ⁻¹	C=C	1650 cm ⁻¹
C≡N	2250 cm ⁻¹	C—O	1100 cm ⁻¹

- The C=O group gives rise to a strong absorption in the region 1820–1600 cm^{-1} . The peak is often the strongest in the spectrum and of medium width. You can't miss it.
- If C=O is present, check the following types. (If it is absent, go to item 3.)

Acids	Is O—H also present? Broad absorption near 3300–2500 cm^{-1} (usually overlaps C—H).
Amides	Is N—H also present? Medium absorption near 3500 cm^{-1} , sometimes a double peak, equivalent halves.
Esters	Is C—O also present? Medium intensity absorptions near 1300–1000 cm^{-1} .
Anhydrides	Have <i>two</i> C=O absorptions near 1810 and 1760 cm^{-1} .
Aldehydes	Is aldehyde C—H present? Two weak absorptions near 2850 cm^{-1} and 2750 cm^{-1} on the right side of C—H absorptions.
Ketones	The preceding five choices have been eliminated.
 - If C=O is absent

Alcohols or Phenols	Check for O—H. <i>Broad</i> absorption near 3600–3300 cm^{-1} . Confirm this by finding C—O near 1300–1000 cm^{-1} .
Amines	Check for N—H. Medium absorption(s) near 3500 cm^{-1} .
Ethers	Check for C—O (and absence of O—H) near 1300–1000 cm^{-1} .
 - Double bonds or aromatic rings or both

C=C is a *weak* absorption near 1650 cm^{-1} .
Medium to strong absorptions in the region 1650–1450 cm^{-1} often imply an aromatic ring.
Confirm the above by consulting the C—H region.
Aromatic and vinyl C—H occur to the left of 3000 cm^{-1} (aliphatic C—H occurs to the right of this value).
 - Triple bonds

C≡N is a medium, sharp absorption near 2250 cm^{-1} .
C≡C is a weak but sharp absorption near 2150 cm^{-1} .
Check also for acetylenic C—H near 3300 cm^{-1} .
 - Nitro groups

Two strong absorptions near 1600–1500 cm^{-1} and 1390–1300 cm^{-1} .
 - Hydrocarbons

None of the above is found.
Main absorptions are in the C—H region near 3000 cm^{-1} .
Very simple spectrum, only other absorptions are near 1450 cm^{-1} and 1375 cm^{-1} .

The beginning student should resist the idea of trying to assign or interpret *every* peak in the spectrum. You simply will not be able to do this. Concentrate first on learning the principal peaks and recognizing their presence or absence. This is best done by carefully studying the illustrative spectra in the section that follows.

NOTE: In describing the shifts of absorption peaks or their relative positions, we have used the phrases “to the left” and “to the right.” This was done to simplify descriptions of peak positions. The meaning is clear because all spectra are conventionally presented left to right from 4000 – 600 cm^{-1} .

25.14 Survey of the Important Functional Groups

A. Alkanes

The spectrum is usually simple, with a few peaks.

C—H Stretch occurs around 3000 cm^{-1} .

1. In alkanes (except strained ring compounds), absorption always occurs to the right of 3000 cm^{-1} .
2. If a compound has vinylic, aromatic, acetylenic, or cyclopropyl hydrogens, the CH absorption is to the left of 3000 cm^{-1} .

CH₂ Methylene groups have a characteristic absorption at approximately 1450 cm^{-1} .

CH₃ Methyl groups have a characteristic absorption at approximately 1375 cm^{-1} .

C—C Stretch—not interpretatively useful—has many peaks.

The spectrum of decane is shown in Figure 25.15.

B. Alkenes

=C—H Stretch occurs to the left of 3000 cm^{-1} .

=C—H Out-of-plane (oop) bending occurs at $1000\text{--}650\text{ cm}^{-1}$.

The C—H oop absorptions often allow you to determine the type of substitution pattern on the double bond, according to the number of absorptions and their positions. The correlation chart in Figure 25.16 shows the positions of these bands.

C=C Stretch $1675\text{--}1600\text{ cm}^{-1}$, often weak.

Conjugation moves C=C stretch to the right.

Symmetrically substituted bonds, as in 2,3-dimethyl-2-butene, do not absorb in the infrared region (no dipole change). Highly substituted double bonds are often vanishingly weak in absorption.

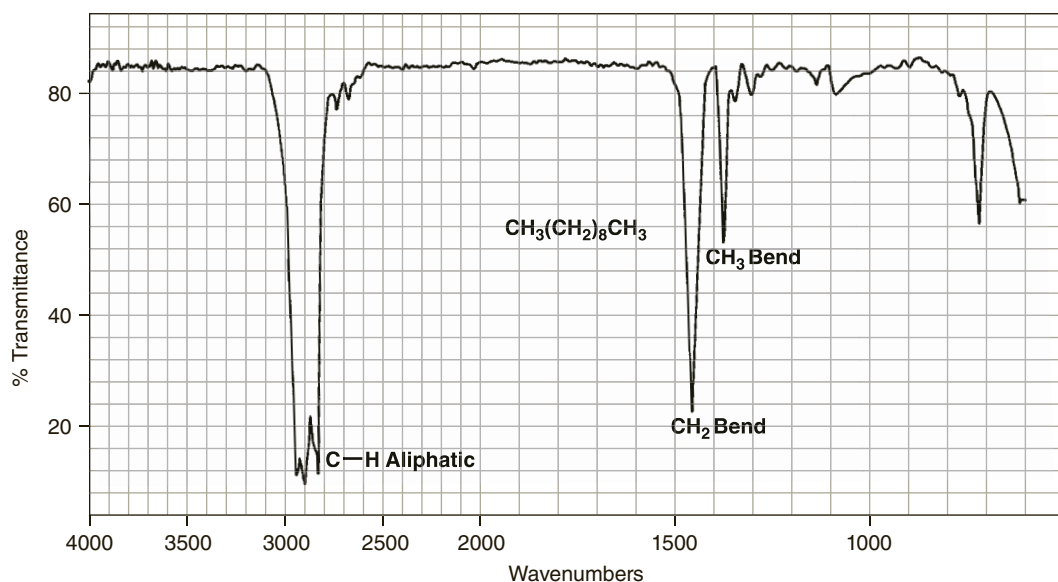


Figure 25.15

Infrared spectrum of decane (neat liquid, salt plates).

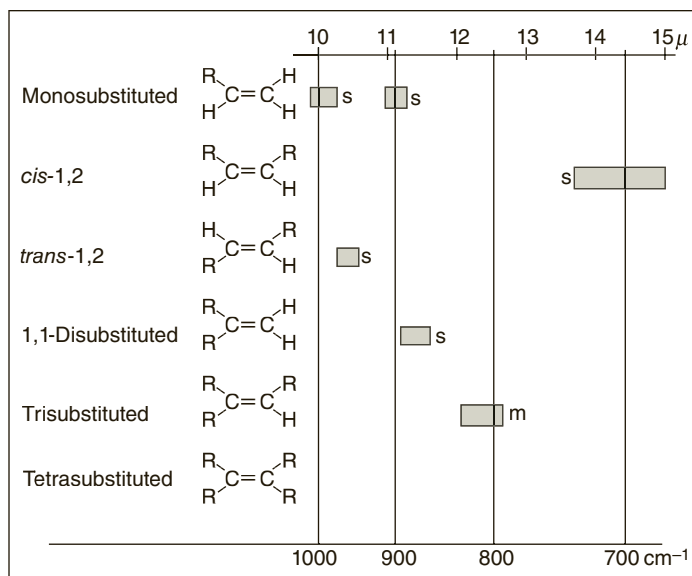


Figure 25.16
The C—H out-of-plane bending vibrations for substituted alkenes.

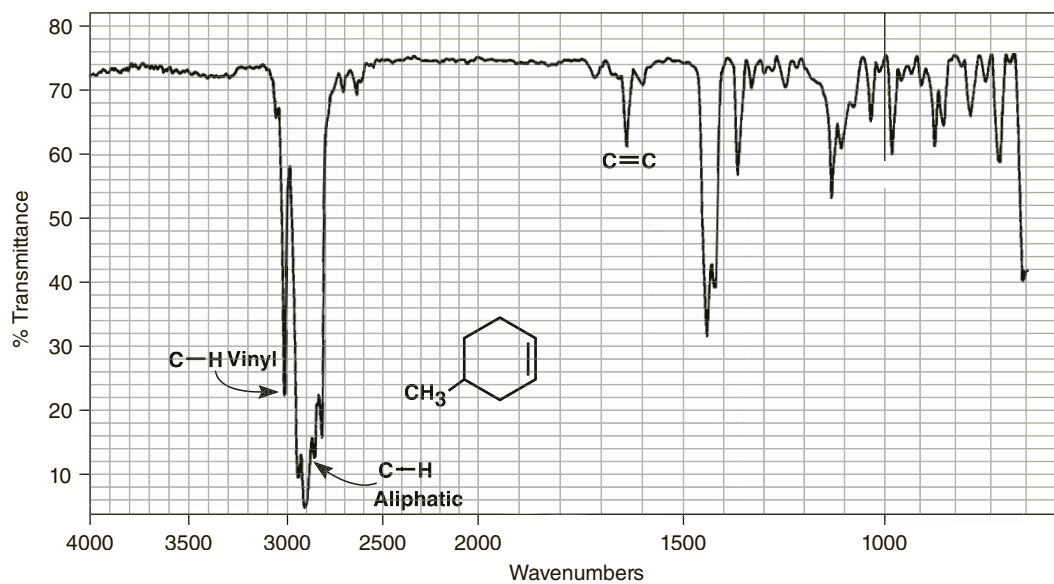


Figure 25.17
Infrared spectrum of 4-methylcyclohexene (neat liquid, salt plates).

The spectra of 4-methylcyclohexene and styrene are shown in Figures 25.17 and 25.18.

C. Aromatic Rings

=C—H Stretch is always to the left of 3000 cm^{-1} .

=C—H Out-of-plane oop bending occurs at $900\text{ to }690\text{ cm}^{-1}$.

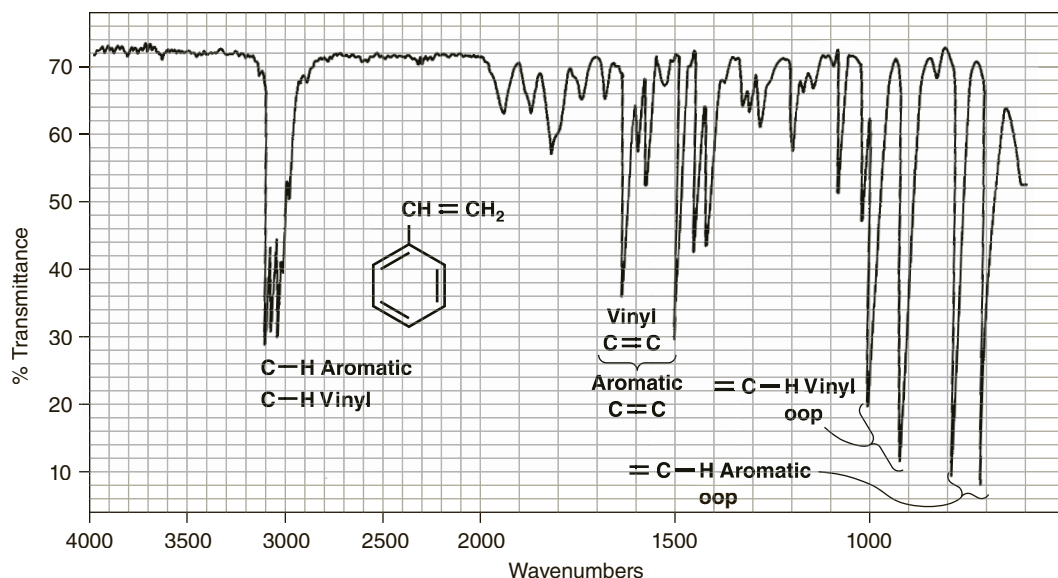


Figure 25.18
Infrared spectrum of styrene (neat liquid, salt plates).

The C—H oop absorptions often allow you to determine the type of ring substitution by their numbers, intensities, and positions. The correlation chart in Figure 25.19A indicates the positions of these bands.

The patterns are generally reliable—they are most reliable for rings with alkyl substituents and least reliable for polar substituents.

Ring Absorptions (C=C). There are often four sharp absorptions that occur in pairs at 1600 cm^{-1} and 1450 cm^{-1} and are characteristic of an aromatic ring. See, for example, the spectra of anisole (Figure 25.23), benzonitrile (Figure 25.26), and methyl benzoate (Figure 25.35).

There are many weak combination and overtone absorptions that appear between 2000 cm^{-1} and 1667 cm^{-1} . The relative shapes and numbers of these peaks can be used to determine whether an aromatic ring is monosubstituted or di-, tri-, tetra-, penta-, or hexa-substituted. Positional isomers can also be distinguished. Because the absorptions are weak, these bands are best observed by using neat liquids or concentrated solutions. If the compound has a high-frequency carbonyl group, this absorption overlaps the weak overtone bands, so no useful information can be obtained from analyzing this region. The various patterns that are obtained in this region are shown in Figure 25.19B.

The spectra of styrene and *o*-dichlorobenzene are shown in Figures 25.18 and 25.20.


D. Alkynes

$\equiv\text{C—H}$ Stretch is usually near 3300 cm^{-1} , sharp peak.


$\text{C}\equiv\text{C}$ Stretch is near 2150 cm^{-1} , sharp peak.

Conjugation moves $\text{C}\equiv\text{C}$ stretch to the right.

Disubstituted or symmetrically substituted triple bonds give either no absorption or weak absorption.

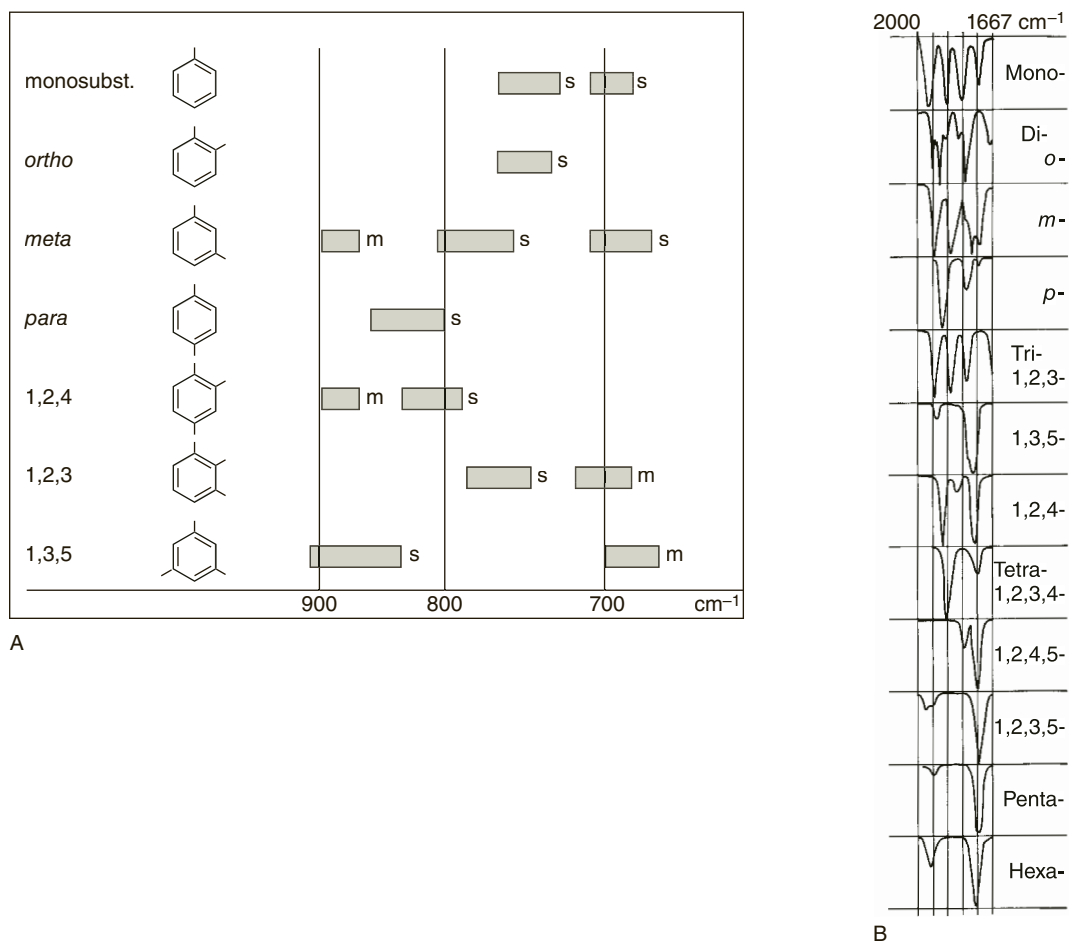**Figure 25.19**

(A) The C—H out-of-plane bending vibrations for substituted benzenoid compounds. (B) The 2000–1667 cm^{-1} region for substituted benzenoid compounds. (From Dyer, J.R. *Applications of Absorption Spectroscopy of Organic Compounds*, Prentice Hall: Englewood Cliffs, NJ, 1965.)

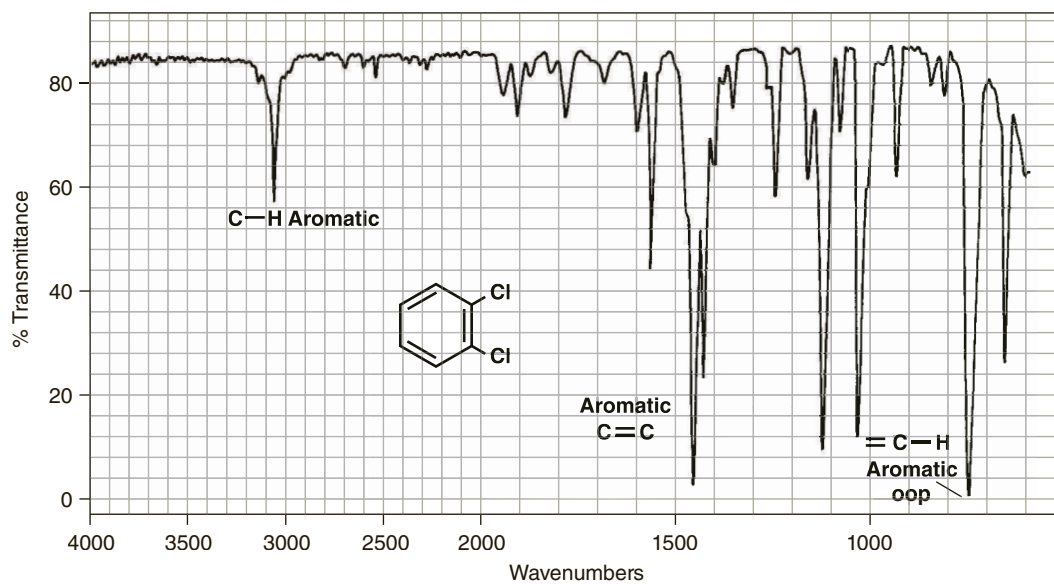**Figure 25.20**

Infrared spectrum of *o*-dichlorobenzene (neat liquid, salt plates).

E. Alcohols and Phenols

- O—H Stretch is a sharp peak at $3650\text{--}3600\text{ cm}^{-1}$ if no hydrogen bonding takes place. (This is usually observed only in dilute solutions.)
If there is hydrogen bonding (usual in neat or concentrated solutions), the absorption is *broad* and occurs more to the right at $3500\text{--}3200\text{ cm}^{-1}$, sometimes overlapping C—H stretch absorptions.
- C—O Stretch is usually in the range of $1300\text{--}1000\text{ cm}^{-1}$.
Phenols are like alcohols. The 2-naphthol shown in Figure 25.21 has some molecules hydrogen bonded and some free. The spectrum of 4-methylcyclohexanol is shown in Figure 25.22. This alcohol, which was determined neat, would also have had a free OH spike to the left of this hydrogen-bonded band if it had been determined in dilute solution.

F. Ethers

- C—O The most prominent band is due to C—O stretch at $1300\text{--}1000\text{ cm}^{-1}$.
Absence of C=O and O—H bands is required to be sure the C—O stretch is not due to an alcohol or ester. Phenyl and vinyl ethers are found in the left portion of the range, aliphatic ethers in the right. (Conjugation with the oxygen moves the absorption to the left.)

The spectrum of anisole is shown in Figure 25.23.

G. Amines

- N—H Stretch occurs in the range of $3500\text{--}3300\text{ cm}^{-1}$.
Primary amines have *two* bands typically 30 cm^{-1} apart.
Secondary amines have one band, often vanishingly weak.
Tertiary amines have no NH stretch.

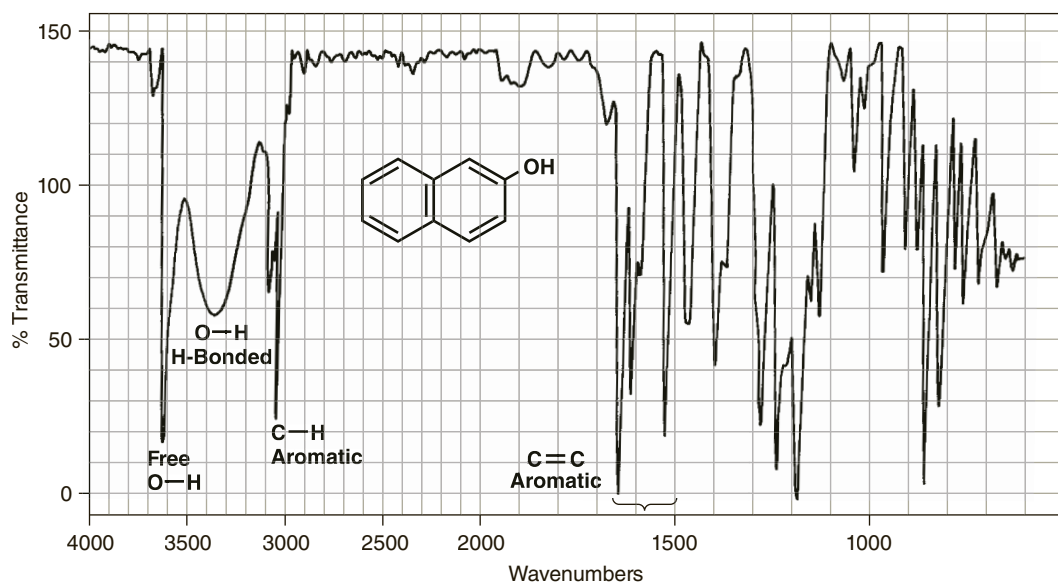


Figure 25.21
Infrared spectrum of 2-naphthol showing both free and hydrogen-bonded OH (CHCl_3 solution).

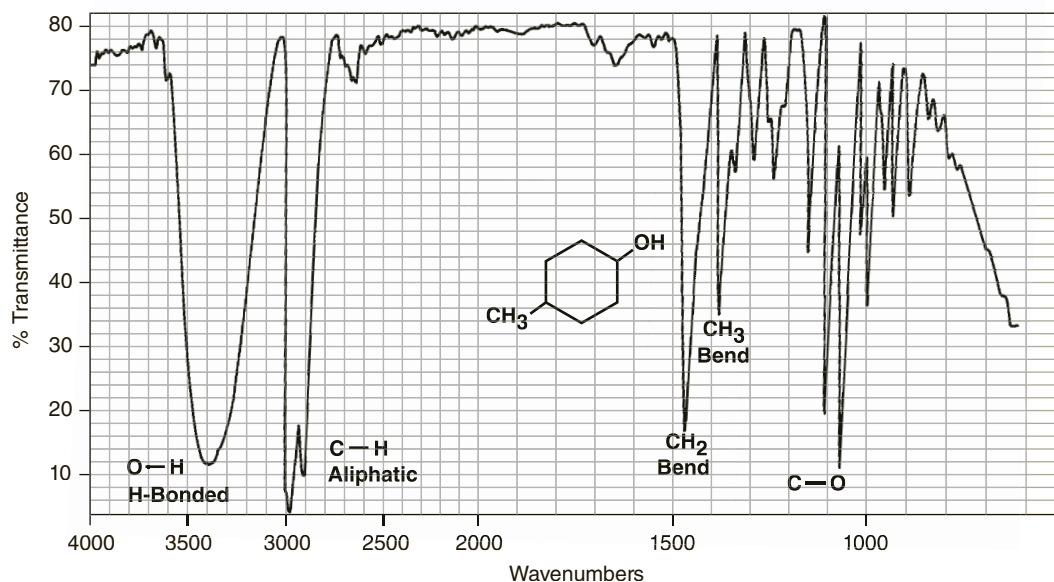


Figure 25.22
Infrared spectrum of 4-methylcyclohexanol (neat liquid, salt plates).

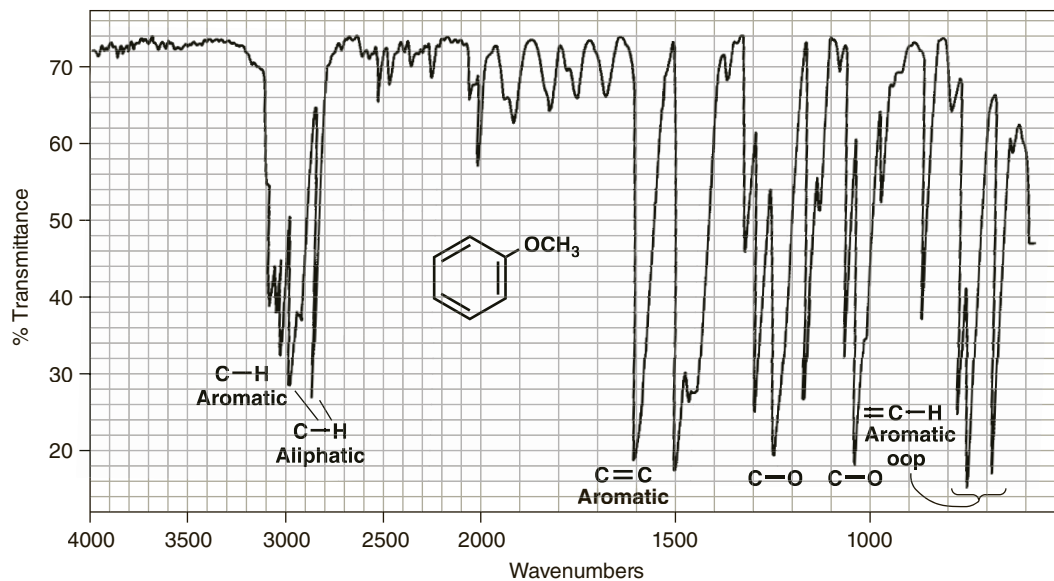


Figure 25.23
Infrared spectrum of anisole (neat liquid, salt plates).

C—N Stretch is weak and occurs in the range of 1350–1000 cm^{-1} .

N—H Scissoring bending mode occurs in the range of 1640–1560 cm^{-1} (broad).

An oop bending absorption can sometimes be observed at about 800 cm^{-1} .

The spectrum of *n*-butylamine is shown in Figure 25.24.

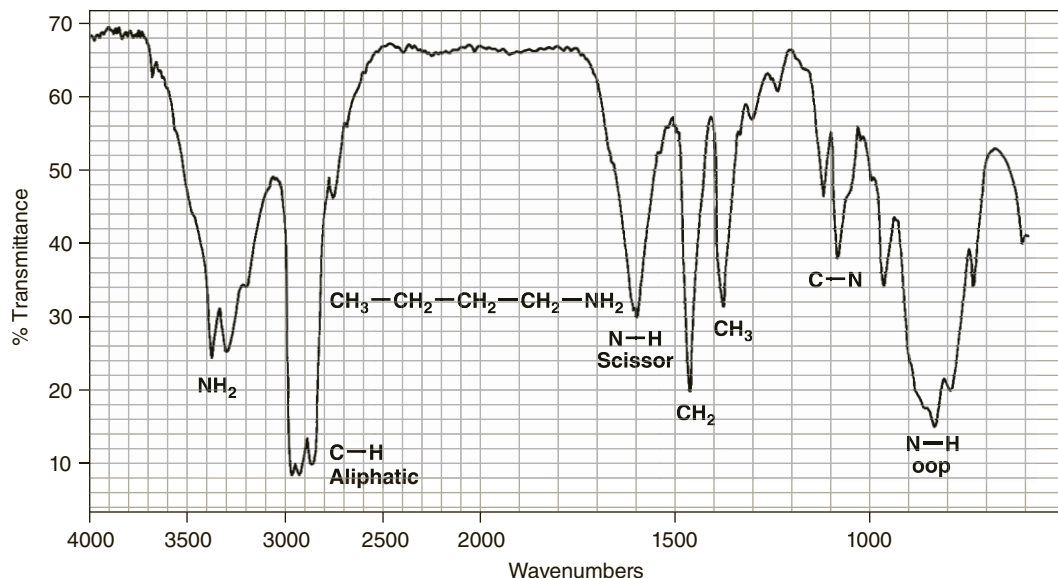


Figure 25.24
Infrared spectrum of *n*-butylamine (neat liquid, salt plates).

H. Nitro Compounds

$\text{N}=\text{O}$ Stretch is usually two strong bands at $1600\text{--}1500\text{ cm}^{-1}$ and $1390\text{--}1300\text{ cm}^{-1}$.

The spectrum of nitrobenzene is shown in Figure 25.25.

I. Nitriles

$\text{C}\equiv\text{N}$ Stretch is a sharp absorption near 2250 cm^{-1} .
Conjugation with double bonds or aromatic rings moves the absorption to the right.

The spectrum of benzonitrile is shown in Figure 25.26.

J. Carbonyl Compounds

The carbonyl group is one of the most strongly absorbing groups in the infrared region of the spectrum. This is mainly due to its large dipole moment. It absorbs in a variety of compounds (aldehydes, ketones, acids, esters, amides, anhydrides, and acid chlorides) in the range of $1850\text{--}1650\text{ cm}^{-1}$. In Figure 25.27, the normal values for the various types of carbonyl groups are compared. In the sections that follow, each type is examined separately.

K. Aldehydes

$\text{C}=\text{O}$ Stretch at approximately 1725 cm^{-1} is normal.
Aldehydes *seldom* absorb to the left of this value.
Conjugation moves the absorption to the right.

$\text{C}-\text{H}$ Stretch, aldehyde hydrogen ($-\text{CHO}$), consists of *weak* bands at about 2750 cm^{-1} and 2850 cm^{-1} . Note that the CH stretch in alkyl chains does not usually extend this far to the right.

The spectrum of an unconjugated aldehyde, nonanal, is shown in Figure 25.28, and the conjugated aldehyde, benzaldehyde, is shown in Figure 25.29.

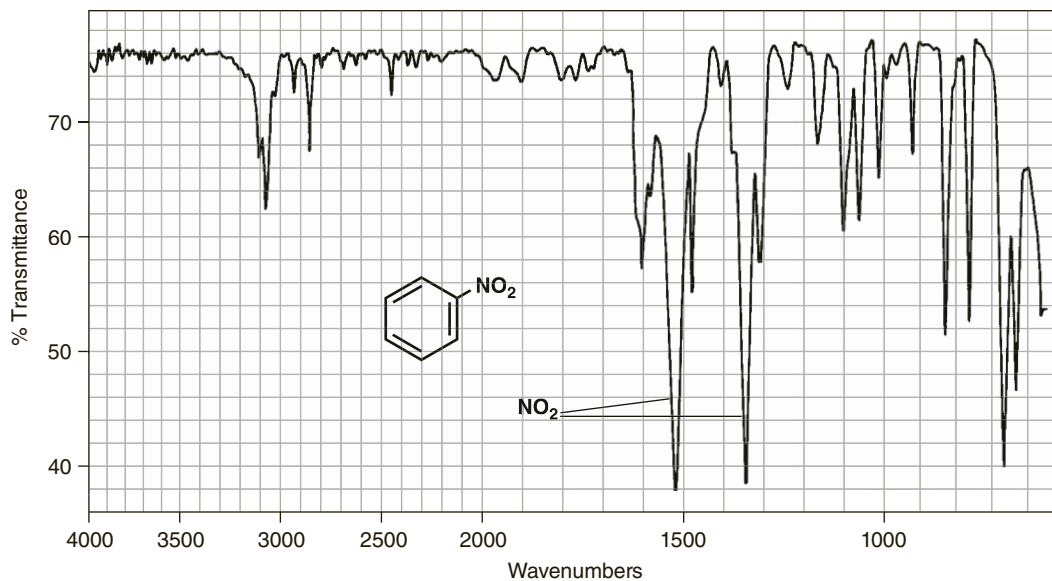


Figure 25.25
Infrared spectrum of nitrobenzene (neat liquid, salt plates).

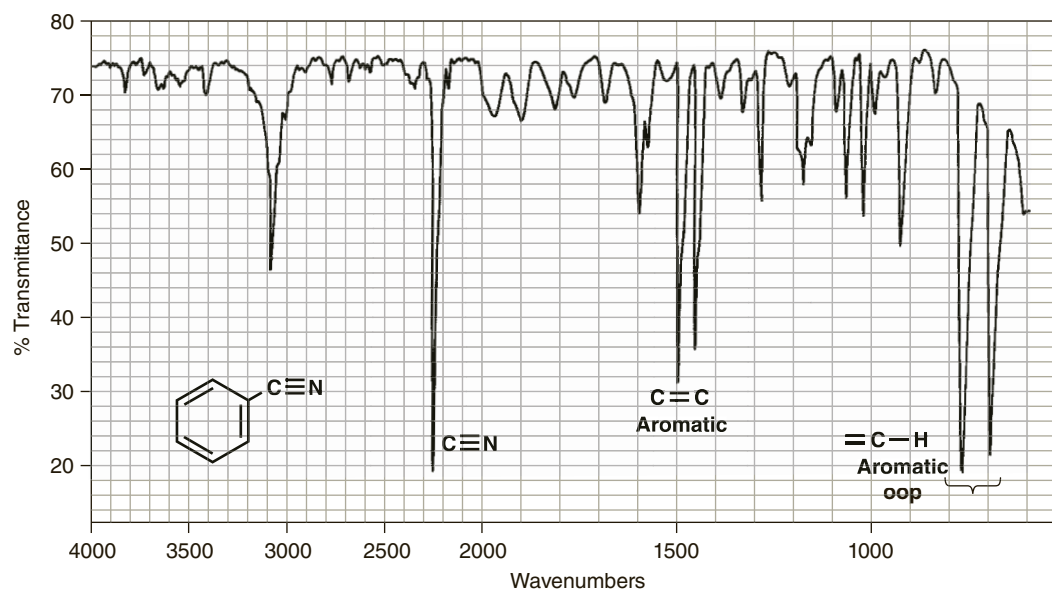


Figure 25.26
Infrared spectrum of benzonitrile (neat liquid, salt plates).

←————— cm ⁻¹ —————→							
1810	1800	1760	1735	1725	1715	1710	1675
Anhydride (band 1)	Acid chloride	Anhydride (band 2)	Ester	Aldehyde	Ketone	Carboxylic acid	Amide

Figure 25.27
Normal base values for the C=O stretching vibrations for carbonyl groups.

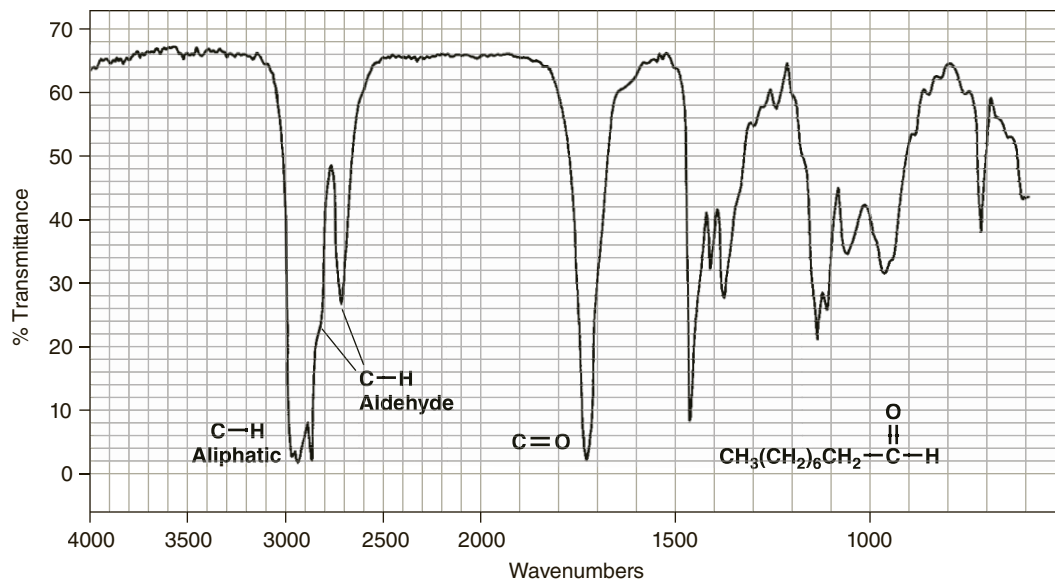


Figure 25.28
Infrared spectrum of nonanal (neat liquid, salt plates).

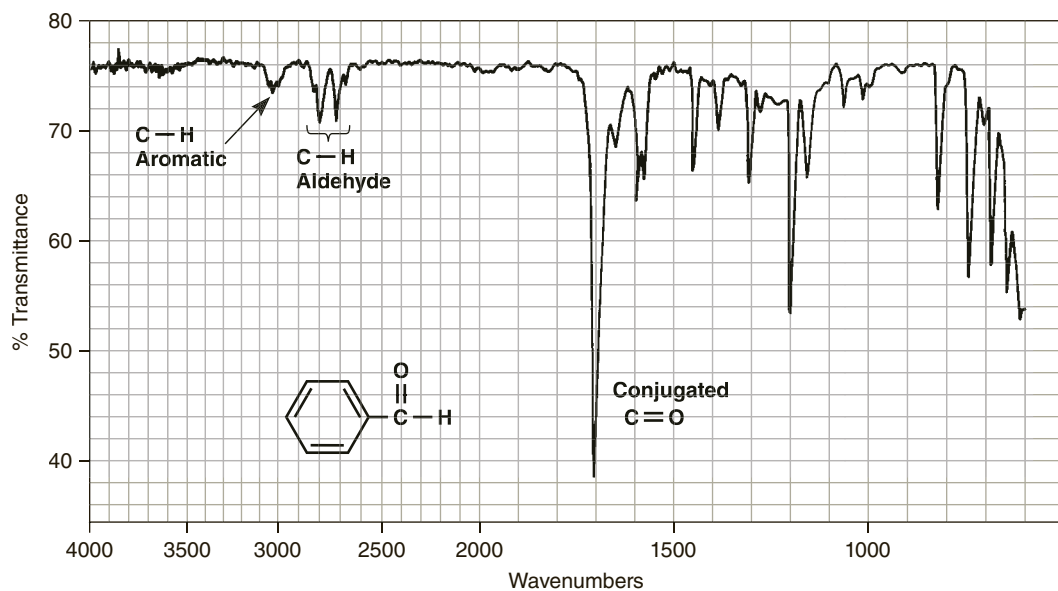


Figure 25.29
Infrared spectrum of benzaldehyde (neat liquid, salt plates).

L. Ketones

- $\text{C}=\text{O}$ Stretch at approximately at 1715 cm^{-1} is normal.
 Conjugation moves the absorption to the right.
 Ring strain moves the absorption to the left in cyclic ketones
 (see Figure 25.30).

The spectra of methyl isopropyl ketone and mesityl oxide are shown in Figures 25.14 and 25.31. The spectrum of camphor, shown in Figure 25.32, has a carbonyl group that has been shifted to a higher frequency because of ring strain (1745 cm^{-1}).

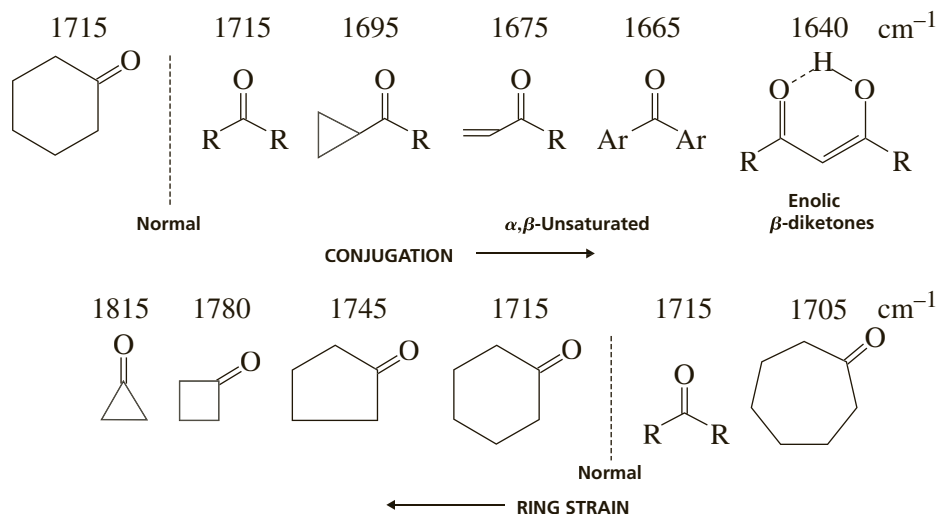


Figure 25.30
Effects of conjugation and ring strain on carbonyl frequencies in ketones.

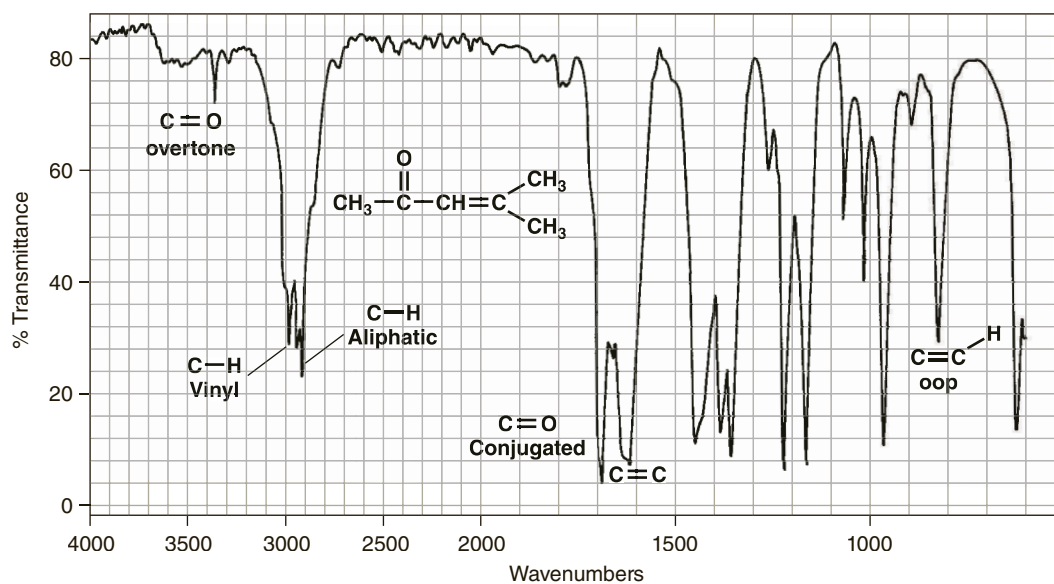


Figure 25.31
Infrared spectrum of mesityl oxide (neat liquid, salt plates).

M. Acids

O-H Stretch, usually *very broad* (strongly hydrogen-bonded) at 3300–2500 cm^{-1} , often interferes with C-H absorptions.

C=O Stretch, broad, 1730–1700 cm^{-1} .
Conjugation moves the absorption to the right.

C-O Stretch, in the range of 1320–1210 cm^{-1} , is strong.

The spectrum of benzoic acid is shown in Figure 25.33.

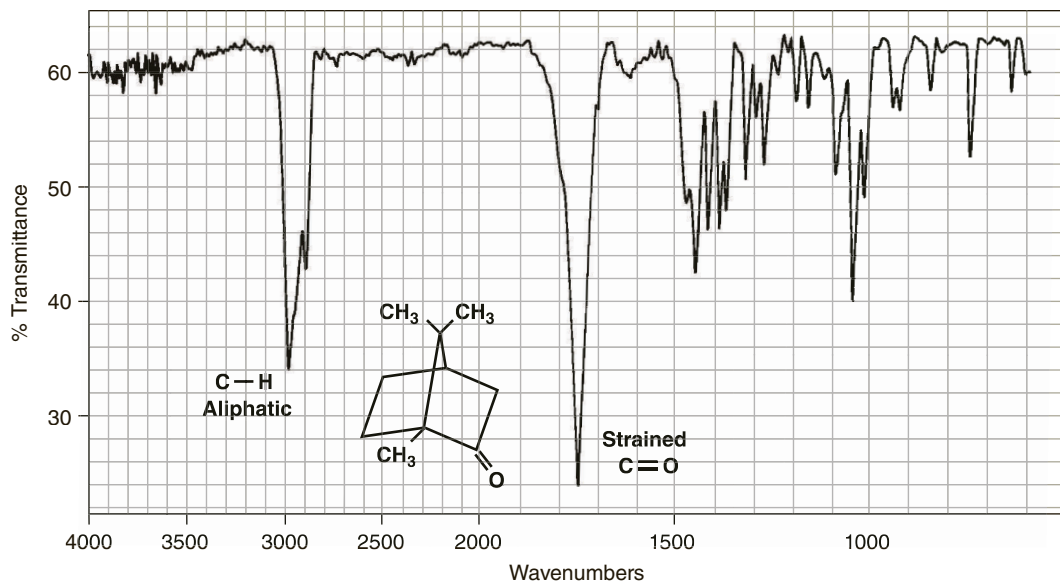


Figure 25.32
Infrared spectrum of camphor (KBr pellet).


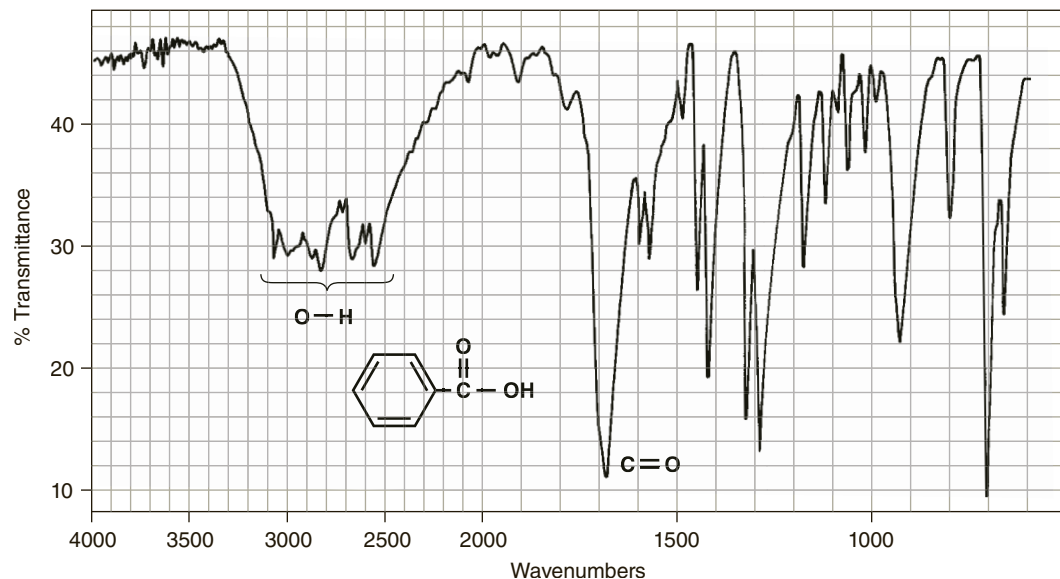


Figure 25.33
Infrared spectrum of benzoic acid (KBr pellet).

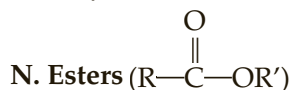

$\text{C}=\text{O}$ Stretch occurs at about 1735 cm^{-1} in normal esters.

1. Conjugation in the R part moves the absorption to the right.
2. Conjugation with the O in the R' part moves the absorption to the left.
3. Ring strain (lactones) moves the absorption to the left.

$\text{C}-\text{O}$ Stretch, two bands or more, one stronger than the others, is in the range of $1300\text{--}1000\text{ cm}^{-1}$.

The spectrum of an unconjugated ester, isopentyl acetate, is shown in Figure 25.34 ($\text{C}=\text{O}$ appears at 1740 cm^{-1}). A conjugated ester, methyl benzoate, is shown in Figure 25.35 ($\text{C}=\text{O}$ appears at 1720 cm^{-1}).

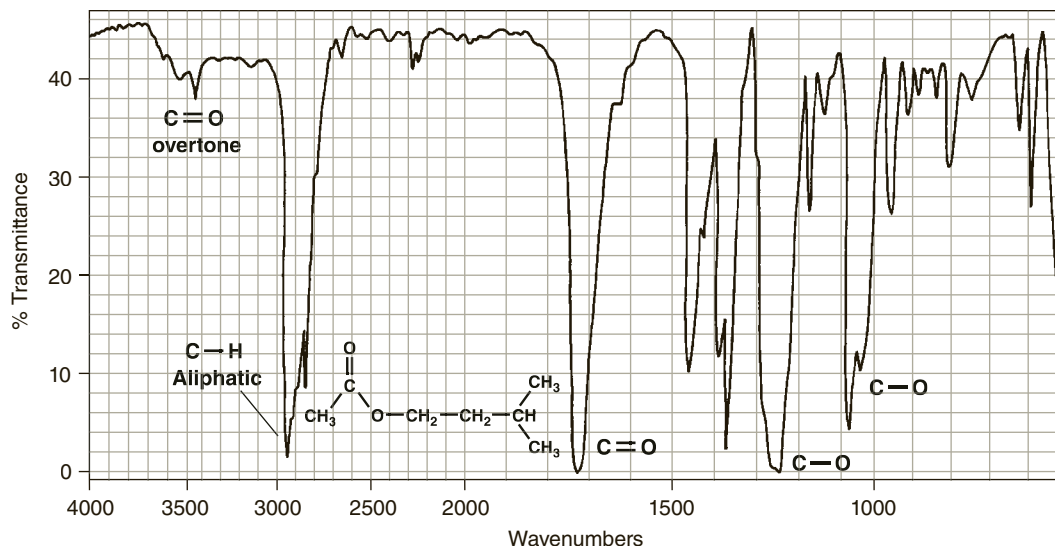


Figure 25.34
Infrared spectrum of isopentyl acetate (neat liquid, salt plates).

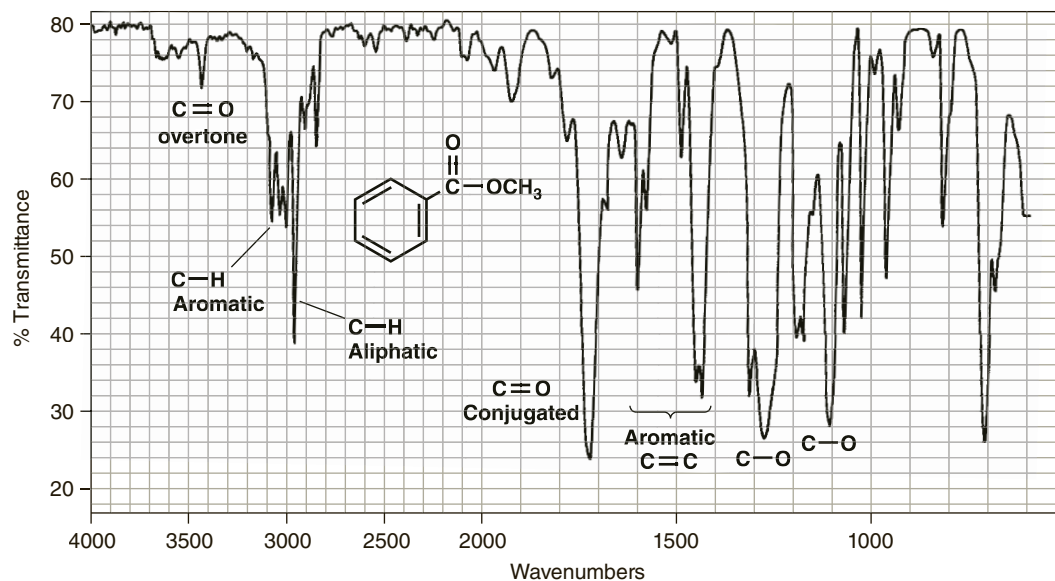


Figure 25.35
Infrared spectrum of methyl benzoate (neat liquid, salt plates).

O. Amides

C=O Stretch is at approximately $1700\text{--}1640\text{ cm}^{-1}$.

Conjugation and ring size (lactams) have the usual effects.

N—H Stretch (if monosubstituted or unsubstituted) is at $3500\text{--}3100\text{ cm}^{-1}$.

Unsubstituted amides have two bands (—NH_2) in this region.

N—H Bending around $1640\text{--}1550\text{ cm}^{-1}$.

The spectrum of benzamide is shown in Figure 25.36.

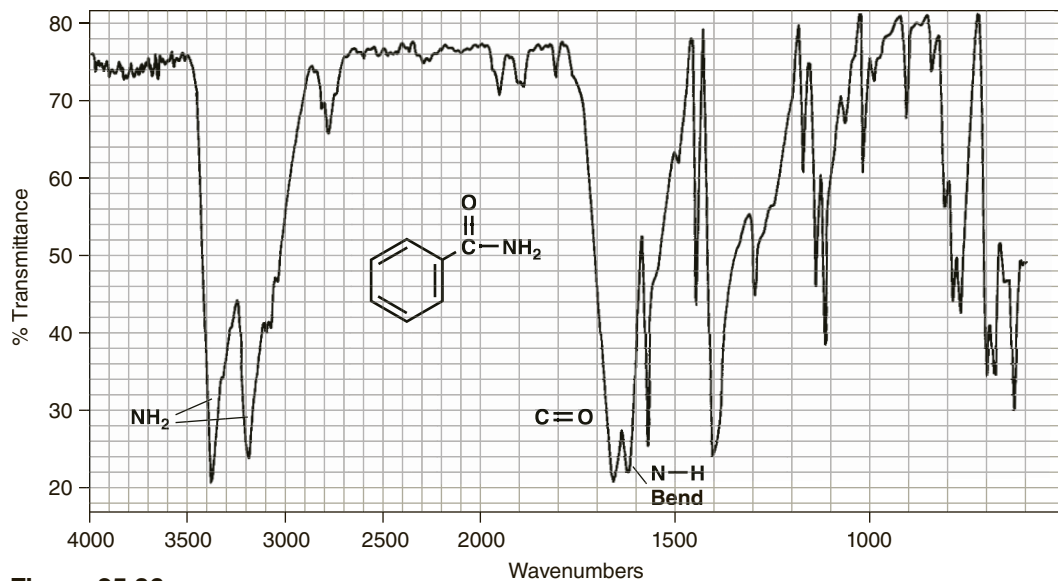


Figure 25.36
Infrared spectrum of benzamide (solid phase, KBr).

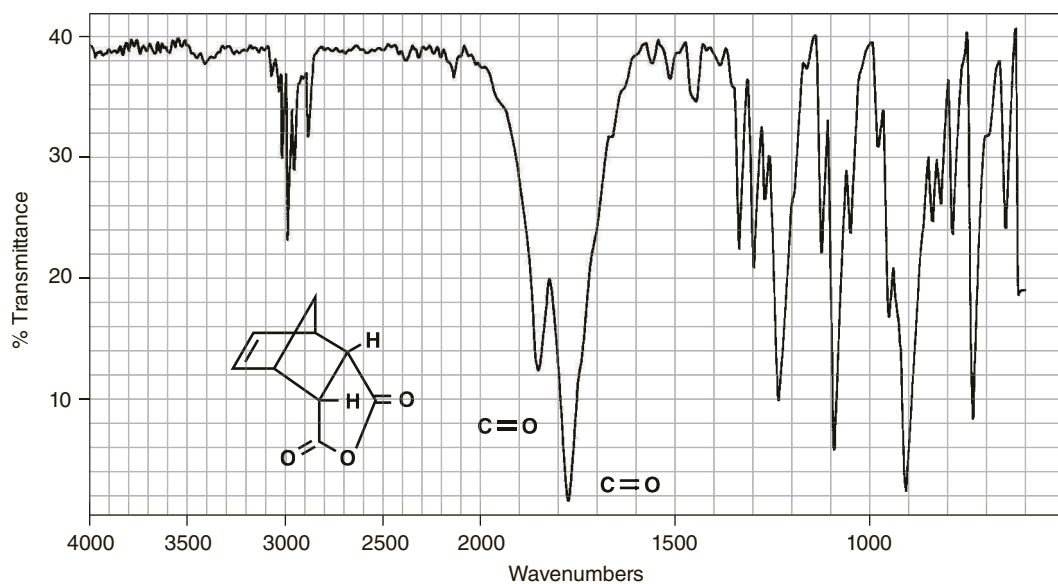


Figure 25.37
Infrared spectrum of *cis*-norbornene-5,6-*endo*-dicarboxylic anhydride (KBr pellet).

P. Anhydrides

- $\text{C}=\text{O}$ Stretch always has *two* bands: $1830\text{--}1800\text{ cm}^{-1}$ and $1775\text{--}1740\text{ cm}^{-1}$.
 Unsaturation moves the absorptions to the right.
 Ring strain (cyclic anhydrides) moves the absorptions to the left.
- $\text{C}-\text{O}$ Stretch is at $1300\text{--}900\text{ cm}^{-1}$. The spectrum of *cis*-norbornene-5,6-*endo*-dicarboxylic anhydride is shown in Figure 25.37.

Q. Acid Chlorides

- C=O Stretch occurs in the range 1810–1775 cm^{-1} in unconjugated chlorides. Conjugation lowers the frequency to 1780–1760 cm^{-1} .
- C—O Stretch occurs in the range 730–550 cm^{-1} .

R. Halides

It is often difficult to determine either the presence or the absence of a halide in a compound by infrared spectroscopy. The absorption bands cannot be relied on, especially if the spectrum is being determined with the compound dissolved in CCl_4 or CHCl_3 solution.

- C—F Stretch, 1350–960 cm^{-1} .
- C—Cl Stretch, 850–500 cm^{-1} .
- C—Br Stretch, to the right of 667 cm^{-1} .
- C—I Stretch, to the right of 667 cm^{-1} .

The spectra of the solvents, carbon tetrachloride and chloroform, are shown in Figures 25.8 and 25.9, respectively.













REFERENCES

- Bellamy, L. J. *The Infra-Red Spectra of Complex Molecules*, 3rd ed.; Methuen: New York, 1975.
- Colthup, N. B.; Daly, L. H.; Wiberly, S. E. *Introduction to Infrared and Raman Spectroscopy*, 3rd ed. Academic Press: San Diego, CA, 1990.
- Dyer, J. R. *Applications of Absorption Spectroscopy of Organic Compounds*; Prentice-Hall: Englewood Cliffs, NJ, 1965.
- Lin-Vien, D.; Colthup, N. B.; Fateley, W. G.; Grasselli, J. G. *Infrared and Raman Characteristic Frequencies of Organic Molecules*; Academic Press: San Diego, CA, 1991.
- Nakanishi, K.; Solomon, P. H. *Infrared Absorption Spectroscopy*, 2nd ed.; Holden-Day: San Francisco, 1977.
- Pavia, D. L.; Lampman, G. M.; Kriz, G. S.; Vyryan, J. R. *Introduction to Spectroscopy: A Guide for Students of Organic Chemistry*, 5th ed. Cengage Learning: Stamford, CT, 2015.
- Silverstein, R. M.; Webster, F. X.; Kiemle, D. J.; Bryce, D. L. *Spectrometric Identification of Organic Compounds*, 8th ed.; Wiley & Sons: New York, 2015.

PROBLEMS

- Comment on the suitability of running the infrared spectrum under each of the following conditions. If there is a problem with the conditions given, provide a suitable alternative method.
 - A neat spectrum of liquid with a boiling point of 150°C is determined using salt plates.
 - A neat spectrum of a liquid with a boiling point of 35°C is determined using salt plates.
 - A KBr pellet is prepared with a compound that melts at 200°C.
 - A KBr pellet is prepared with a compound that melts at 30°C.
 - A solid aliphatic hydrocarbon compound is determined as a Nujol mull.

- f. Silver chloride plates are used to determine the spectrum of aniline.
 g. Sodium chloride plates are selected to run the spectrum of a compound that contains some water.
2. Indicate how you could distinguish between the following pairs of compounds by using infrared spectroscopy.

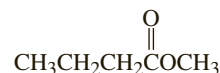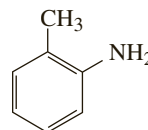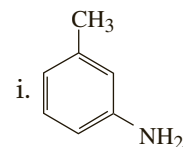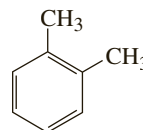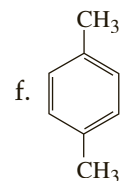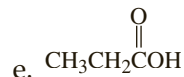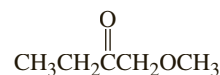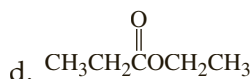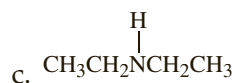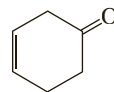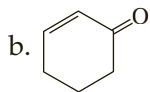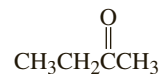

Nuclear Magnetic Resonance Spectroscopy (Proton NMR)

Nuclear magnetic resonance (NMR) spectroscopy is an instrumental technique that allows the number, type, and relative positions of certain atoms in a molecule to be determined. This type of spectroscopy applies only to those atoms that have nuclear magnetic moments because of their nuclear spin properties. Although many atoms meet this requirement, hydrogen atoms (^1H) are of the greatest interest to the organic chemist. Atoms of the ordinary isotopes of carbon (^{12}C) and oxygen (^{16}O) do not have nuclear magnetic moments, and ordinary nitrogen atoms (^{14}N), although they do have magnetic moments, generally fail to show typical NMR behavior for other reasons. The same is true of the halogen atoms, except for fluorine (^{19}F), which does show active NMR behavior. Of the atoms mentioned here, the hydrogen nucleus (^1H) and carbon-13 nucleus (^{13}C) are the most important to organic chemists. Proton (^1H) NMR is discussed here and carbon (^{13}C) NMR is described in Technique 27.

Nuclei of NMR-active atoms placed in a magnetic field can be thought of as tiny bar magnets. In hydrogen, which has two allowed nuclear spin states ($+\frac{1}{2}$ and $-\frac{1}{2}$), either the nuclear magnets of individual atoms can be aligned with the magnetic field (spin $+\frac{1}{2}$), or they can be opposed to it (spin $-\frac{1}{2}$). A slight majority of the nuclei are aligned with the field, because this spin orientation constitutes a slightly lower-energy spin state. If radio-frequency waves of the appropriate energy are supplied, nuclei aligned with the field can absorb this radiation and reverse their direction of spin or become reoriented so that the nuclear magnet opposes the applied magnetic field (see Figure 26.1).

The frequency of radiation required to induce spin conversion is a direct function of the strength of the applied magnetic field. When a spinning hydrogen nucleus is placed in a magnetic field, the nucleus begins to precess with angular frequency ω , much like a child's toy top. This precessional motion is depicted in Figure 26.2. The angular frequency of nuclear precession ω increases as the strength of the applied magnetic field is increased. The radiation that must be supplied to induce spin conversion in a hydrogen nucleus of spin $+\frac{1}{2}$ must have a frequency that just matches the angular precessional frequency ω . This is called the resonance condition, and spin conversion is said to be a resonance process.

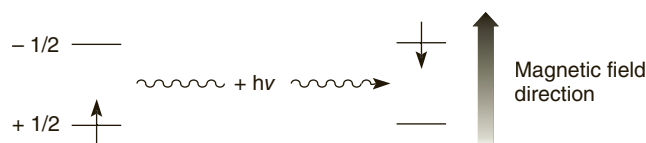


Figure 26.1
The NMR absorption process.

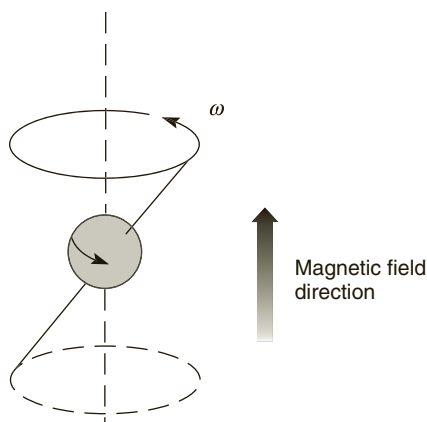


Figure 26.2
Precessional motion of a spinning nucleus in an applied magnetic field.

For the average proton (hydrogen atom), if a magnetic field of approximately 1.4 tesla is applied, radio-frequency radiation of 60 MHz is required to induce a spin transition.¹ Fortunately, the magnetic field strength required to induce the various protons in a molecule to absorb 60-MHz radiation varies from proton to proton within the molecule and is a sensitive function of the immediate *electronic* environment of each proton. The proton nuclear magnetic resonance spectrometer supplies a basic radio-frequency radiation of 60 MHz to the sample being measured and *increases* the strength of the applied magnetic field over a range of several parts per million from the basic field strength. As the field increases, various protons come into resonance (absorb 60-MHz energy), and a resonance signal is generated for each proton. An NMR spectrum is a plot of the strength of the magnetic field versus the intensity of the absorptions. A typical 60-MHz NMR spectrum is shown in Figure 26.3.

Modern FT-NMR instruments produce the same type of NMR spectrum just described, even though they do it by a different method. See your lecture textbook for a discussion of the differences between classic CW instruments and modern FT-NMR instruments. Fourier transform spectrometers operating at magnetic field strengths of at least 7.1 tesla and at spectrometer frequencies of 300 MHz and above allow chemists to obtain both the proton and carbon NMR spectra on the same sample.

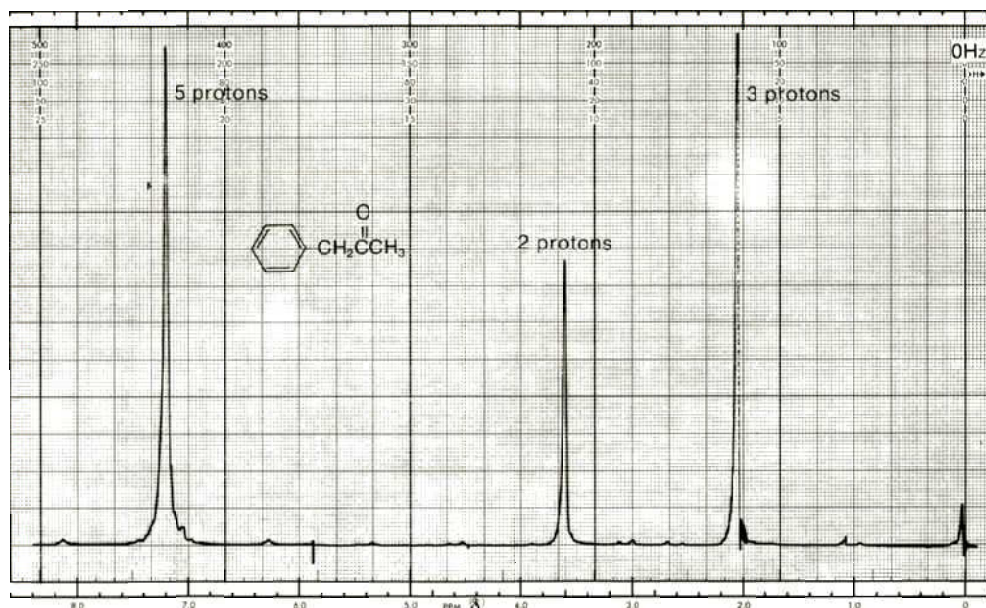


Figure 26.3
Nuclear magnetic resonance spectrum of phenylacetone (the absorption peak at the far right is caused by the added reference substance tetramethylsilane).

¹ Most modern instruments (FT-NMR instruments) use higher fields than described here and operate differently. The classical 60-MHz continuous wave (CW) instrument is used here as a simple example.

PART A. PREPARING A SAMPLE FOR NMR SPECTROSCOPY

The NMR sample tubes used in most instruments are approximately 0.5 cm × 18 cm in overall dimension and are fabricated of uniformly thin glass tubing. These tubes are very fragile and expensive, so care must be taken to avoid breaking the tubes.

CAUTION



NMR tubes are made out of very thin glass and break easily. Never place the cap on tightly, and take special care when removing it.

To prepare the solution, you must first choose the appropriate solvent. The solvent should not have NMR absorption peaks of its own; that is, it should contain no protons. Carbon tetrachloride (CCl_4) fits this requirement and can be used in some instruments. However, because FT-NMR spectrometers require deuterium to stabilize (lock) the field, organic chemists usually use deuterated chloroform (CDCl_3) as a solvent. This solvent dissolves most organic compounds and is relatively inexpensive. You can use this solvent with any NMR instrument. You should not use normal chloroform CHCl_3 because the solvent contains a proton. Deuterium ^2H does not absorb in the proton region and is thus “invisible,” or not seen, in the proton NMR spectrum. Use deuterated chloroform to dissolve your sample, unless you are instructed to use another solvent, such as deuterated derivatives of water, acetone, or dimethylsulfoxide.

26.1 Routine Sample Preparation Using Deuterated Chloroform

1. Most organic liquids and low-melting solids will dissolve in deuterated chloroform. However, you should first determine whether your sample will dissolve in ordinary CHCl_3 before using the deuterated solvent. If your sample does not dissolve in chloroform, consult your instructor about a possible alternative solvent, or consult Section 26.2.

CAUTION



Chloroform, deuterated chloroform, and carbon tetrachloride are all toxic solvents. In addition, they may be carcinogenic substances.

2. If you are using an FT-NMR spectrometer, add 30 mg (0.030 g) of your liquid or solid sample to a tared conical vial or test tube. Use a Pasteur pipette to transfer a liquid or a spatula to transfer a solid. Non-FT instruments usually require a more concentrated solution in order to obtain an adequate spectrum. Typically, a 10–30% sample concentration (weight/weight) is used.
3. Transfer about 0.5 mL of the deuterated chloroform with a clean, dry Pasteur pipette to your sample. Swirl the test tube or conical vial to help dissolve the sample. At this point, the sample should have completely dissolved. Add a little more solvent, if necessary, to dissolve the sample fully.
4. Transfer the solution to the NMR tube using a clean, dry Pasteur pipette. Be careful when transferring the solution to avoid breaking the edge of the fragile NMR tube. It is best to hold the NMR tube and the container with the solution in the same hand when making the transfer.
5. Once the solution has been transferred to the NMR tube, use a clean pipette to add enough deuterated chloroform to bring the total solution height to about 50 mm (see Figure 26.4). In some cases, you will need to add a small amount

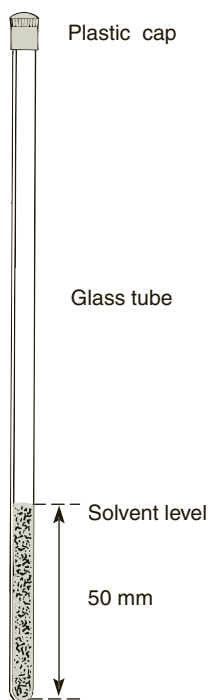


Figure 26.4
An NMR sample tube.

of tetramethylsilane (TMS) as a reference substance (see Section 26.3). Check with your instructor to see if you need to add TMS to your sample. Deuterated chloroform has a small amount of CHCl_3 impurity, which gives rise to a low-intensity peak in the NMR spectrum at 7.27 parts per million (ppm). This impurity may also help you to “reference” your spectrum.

6. Cap the NMR tube. Do this firmly but not too tightly. If you jam the cap on, you may have trouble removing it later without breaking the end off of the very thin glass tube. Make sure that the cap is on straight. Invert the NMR tube several times to mix the contents.
7. You are now ready to record the NMR spectrum of your sample. Insert the NMR tube into its holder and adjust its depth by using the gauge provided to you.

Cleaning the NMR Tube

1. Carefully uncap the tube so that you do not break it. Turn the tube upside down and hold it vertically over a beaker. Shake the tube up and down gently so that its contents empty into the beaker.
2. Partially refill the NMR tube with acetone using a Pasteur pipette. Carefully replace the cap and invert the tube several times to rinse it.
3. Remove the cap and drain the tube as before. Place the open tube upside down in a beaker with a Kimwipe or paper towel placed in the bottom of the beaker. Leave the tube standing in this position for at least one laboratory period so that the acetone completely evaporates. Alternatively, you may place the beaker and NMR tube in an oven for at least 2 hours. If you need to use the NMR tube before the acetone has fully evaporated, attach a piece of pressure tubing to the tube and pull a vacuum with an aspirator. After several minutes, the acetone should have fully evaporated. Because acetone contains protons, you must not use the NMR tube until the acetone has evaporated completely.²
4. Once the acetone is evaporated, place the clean tube and its cap (do not cap the tube) in its storage container and place it in your desk. The storage container will prevent the tube from being crushed.

Health Hazards Associated With NMR Solvents

Carbon tetrachloride, chloroform (and chloroform-d), and benzene (and benzene- d_6) are hazardous solvents. Besides being highly toxic, they are suspected carcinogens. In spite of these health problems, these solvents are commonly used in NMR spectroscopy. Deuterated acetone may be a safer alternative. These solvents are used because they contain no protons and are excellent solvents for most organic compounds. Therefore, you must learn to handle these solvents with great care to minimize the hazard. These solvents should be stored either under a hood or in septum-capped bottles. If the bottles have screw caps, a pipette should be attached to each bottle. A recommended way of attaching the pipette is to store it in a test tube taped to the side of the bottle. Septum-capped bottles can be used only by withdrawing the solvent with a hypodermic syringe that has been designated solely for this use. All samples should be prepared under a hood, and solutions should be disposed of in an appropriately designated waste container that is stored under the hood. Wear rubber or plastic gloves when preparing or discarding samples.


² If you can't wait to be sure all acetone has evaporated, you may rinse the tube once or twice with a *very small* amount of CDCl_3 before using it.

26.2 Nonroutine Sample Preparation

Some compounds do not dissolve readily in CDCl_3 . A commercial solvent called Unisol will often dissolve the difficult cases. Unisol is a mixture of CDCl_3 and DMSO-d_6 . Deuterated acetone may also dissolve more polar substances.

With highly polar substances, you may find that your sample will not dissolve in deuterated chloroform or Unisol. If this is the case, you may be able to dissolve the sample in deuterium oxide D_2O . Spectra determined in D_2O often show a small peak at about 5 ppm because of OH impurity. If the sample compound has acidic hydrogens, they may *exchange* with D_2O , leading to the appearance of an OH peak in the spectrum and the *loss* of the original absorption from the acidic proton, owing to the exchanged hydrogen. In many cases, this will also alter the splitting patterns of a compound.

Many solid carboxylic acids do not dissolve in CDCl_3 or even D_2O . In such cases, add a small piece of sodium metal to about 1 mL of D_2O . The acid is then dissolved in this solution. The resulting basic solution enhances the solubility of the carboxylic acid. In such a case, the hydroxyl proton of the carboxylic acid cannot be observed in the NMR spectrum because it exchanges with the solvent. A large DOH peak is observed, however, due to the exchange and the H_2O impurity in the D_2O solvent.

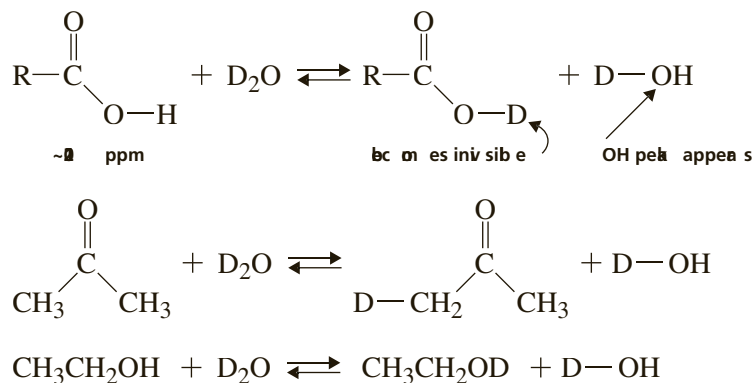

When the above solvents fail, other special solvents can be used. Acetone, acetonitrile, dimethylsulfoxide, pyridine, benzene, and dimethylformamide can be used if you are not interested in the region or regions of the NMR spectrum in which they give rise to absorption. The deuterated (but expensive) analogs of these compounds are also used in special instances (for example, acetone- d_6 , dimethylsulfoxide- d_6 , dimethylformamide- d_7 , and benzene- d_6). If the sample is not sensitive to acid, trifluoroacetic acid (which has no protons with $\delta < 12$) can be used. You must be aware that these solvents often lead to chemical shift values different from those determined in CCl_4 or CDCl_3 . Variations of as much as 0.5–1.0 ppm have been observed. In fact, it is sometimes possible, by switching to pyridine, benzene, acetone, or dimethylsulfoxide as solvents, to separate peaks that overlap when CCl_4 or CDCl_3 solutions are used.


26.3 Reference Substances

To provide the internal reference standard, TMS must be added to the sample solution. This substance has the formula $(\text{CH}_3)_4\text{Si}$. By universal convention, the chemical shifts of the protons in this substance are defined as 0.00 ppm. The spectrum should be shifted so that the TMS signal appears at this position on precalibrated paper.

The concentration of TMS in the sample should range from 1% to 3%. Some people prefer to add 1 to 2 drops of TMS to the sample just before determining

the spectrum. Because TMS has 12 equivalent protons, not much of it needs to be added. A Pasteur pipette or a syringe may be used for the addition. It is far easier to have available in the laboratory a prepared solvent that already contains TMS. Deuterated chloroform and carbon tetrachloride often have TMS added to them. Because TMS is highly volatile (bp 26.5°C), such solutions should be stored, tightly stoppered, in a refrigerator. Tetramethylsilane itself is best stored in a refrigerator as well.

Tetramethylsilane does not dissolve in D₂O. For spectra determined in D₂O, a different internal standard, sodium 2,2-dimethyl-2-silapentane-5-sulfonate, must be used. This standard is water soluble and gives a resonance peak at 0.00 ppm.

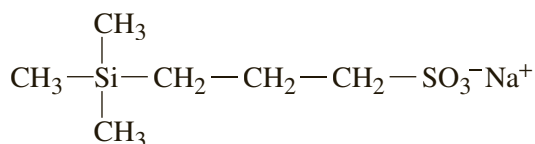

Sodium 2,2-dimethyl-2-silapentane-5-sulfonate (DSS)

PART B. NUCLEAR MAGNETIC RESONANCE (¹H NMR)

26.4 The Chemical Shift

The differences in the applied field strengths at which the various protons in a molecule absorb 60-MHz radiation are extremely small. The different absorption positions amount to a difference of only a few parts per million (ppm) in the magnetic field strength. Because it is experimentally difficult to measure the precise field strength at which each proton absorbs to less than one part in a million, a technique has been developed whereby the *difference* between two absorption positions is measured directly. A standard reference substance is used to achieve this measurement, and the positions of the absorptions of all other protons are measured relative to the values for the reference substance. The reference substance that has been universally accepted is **tetramethylsilane** (CH₃)₄Si, which is also called TMS. The proton resonances in this molecule appear at a higher field strength than the proton resonances in most other molecules, and all the protons of TMS have resonance at the same field strength.

To give the position of absorption of a proton, a quantitative measurement, a parameter called the **chemical shift** (δ), has been defined. One δ unit corresponds to a one-ppm change in the magnetic field strength. To determine the chemical shift value for the various protons in a molecule, the operator determines an NMR spectrum of the molecule with a small quantity of TMS added directly to the sample. That is, both spectra are determined *simultaneously*. The TMS absorption is adjusted to correspond to the δ = 0 ppm position on the recording chart, which is calibrated in δ units, and the δ = 0 values of the absorption peaks for all other protons can be read directly from the chart.

Because the NMR spectrometer increases the magnetic field as the pen moves from left to right on the chart, the TMS absorption appears at the extreme right edge of the spectrum (δ = 0 ppm) or at the *upfield* end of the spectrum. The chart is calibrated in δ units (or ppm), and most other protons absorb at a lower field strength (or *downfield*) from TMS.

The shift from TMS for a given proton depends on the strength of the applied magnetic field. In an applied field of 1.41 tesla, the resonance of a proton is approximately 60 MHz, whereas in an applied field of 2.35 tesla (23,500 gauss), the

resonance appears at approximately 100 MHz. The ratio of the resonance frequencies is the same as the ratio of the two field strengths:

$$\frac{100 \text{ MHz}}{60 \text{ MHz}} = \frac{2.35 \text{ Tesla}}{1.41 \text{ Tesla}} = \frac{23,500 \text{ Gauss}}{14,100 \text{ Gauss}} = \frac{5}{3}$$

Hence, for a given proton, the shift (in hertz) from TMS is five-thirds larger in the 100-MHz range than in the 60-MHz range. This can be confusing for workers trying to compare data if they have spectrometers that differ in the strength of the applied magnetic field. The confusion is easily overcome by defining a new parameter that is independent of field strength—for instance, by dividing the shift in hertz of a given proton by the frequency in megahertz of the spectrometer with which the shift value was obtained. In this manner, a field-independent measure called the **chemical shift** (δ) is obtained:


$$\delta = \frac{\text{(shift in Hz)}}{\text{(spectrometer frequency in MHz)}} \quad (1)$$

The chemical shift in δ units expresses the amount by which a proton resonance is shifted from TMS, in parts per million (ppm), of the spectrometer's basic operating frequency. Values of δ for a given proton are always the same, irrespective of whether the measurement was made at 60 MHz, 100 MHz, or 300 MHz. For instance, at 60 MHz, the shift of the protons in CH_3Br is 162 Hz from TMS; at 100 MHz, the shift is 270 Hz; and at 300 MHz, the shift is 810 Hz. However, all three correspond to the same value of $\delta = 2.70$ ppm:


$$\delta = \frac{162 \text{ Hz}}{60 \text{ MHz}} = \frac{270 \text{ Hz}}{100 \text{ MHz}} = \frac{810 \text{ Hz}}{300 \text{ MHz}} = 2.70 \text{ ppm}$$

26.5 Chemical Equivalence—Integrals

All of the protons in a molecule that are in chemically identical environments often exhibit the same chemical shift. Thus, all of the protons in TMS or all of the protons in benzene, cyclopentane, or acetone have their own respective resonance values all at the same δ value. Each compound gives rise to a single absorption peak in its NMR spectrum. The protons are said to be **chemically equivalent**. On the other hand, molecules that have sets of protons that are chemically distinct from one another may give rise to an absorption peak from each set.

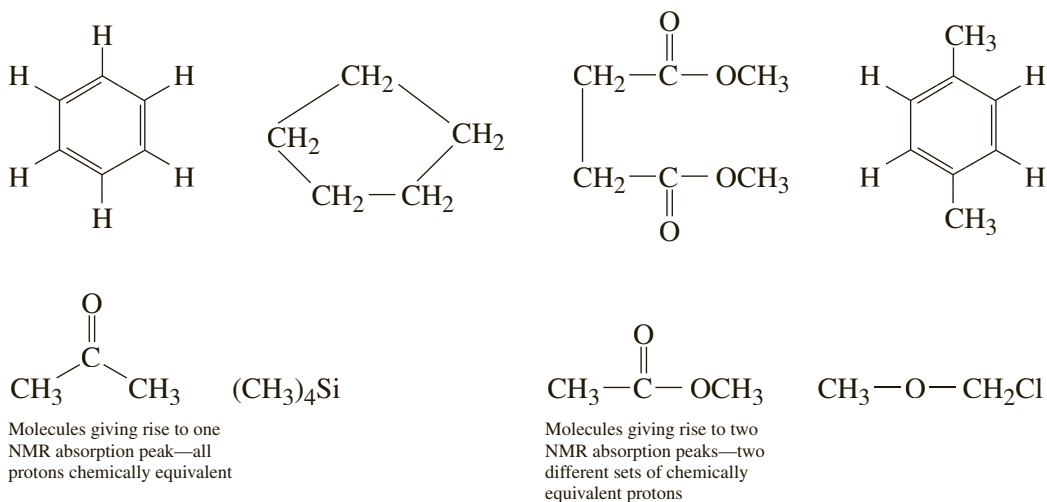

The NMR spectrum given in Figure 26.3 is that of phenylacetone, a compound having *three* chemically distinct types of protons:

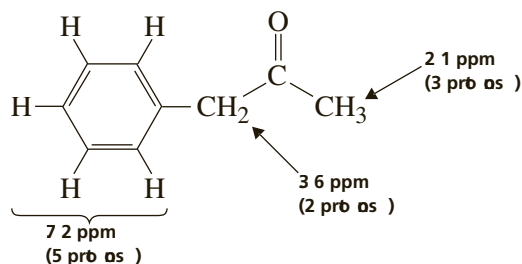

You can immediately see that the NMR spectrum furnishes valuable information on this basis alone. In fact, the NMR spectrum not only can distinguish how many types of protons a molecule has, but also can reveal *how many* of each type are contained within the molecule.

In the NMR spectrum, the area under each peak is proportional to the number of hydrogens generating that peak. Hence, in the case of phenylacetone, the area ratio of the three peaks is 5:2:3, the same as the ratio of the numbers of each type of hydrogen. The NMR spectrometer can electronically “integrate” the area under each peak. It does this by tracing over each peak a vertically rising line, which rises in height by an amount proportional to the area under the peak. Shown in Figure 26.5 is an NMR spectrum of benzyl acetate, with each of the peaks integrated in this way.

It is important to note that the height of the integral line does not give the absolute number of hydrogens; it gives the *relative* numbers of each type of hydrogen. For a given integral to be of any use, there must be a second integral to which it is referred. The benzyl acetate case provides a good example of this. The first integral rises for 55.5 divisions on the chart paper, the second for 22.0 divisions, and the third for 32.5 divisions. These numbers are relative and give the *ratios* of the

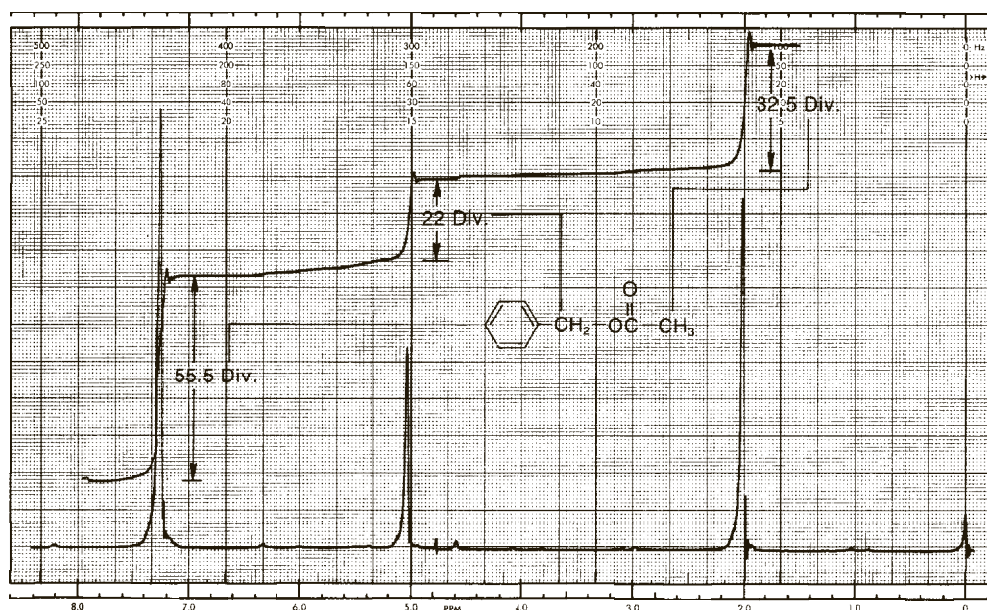


Figure 26.5
Determination of the integral ratios for benzyl acetate.

various types of protons. You can find these ratios by dividing each of the larger numbers by the smallest number:

$$\frac{55.5 \text{ div}}{22.0 \text{ div}} = 2.52 \quad \frac{22.0 \text{ div}}{22.0 \text{ div}} = 1.00 \quad \frac{32.5 \text{ div}}{22.0 \text{ div}} = 1.48$$

Thus, the number ratio of the protons of each type is 2.52:1.00:1.48. If you assume that the peak at 5.1 ppm is really caused by two hydrogens and that the integrals are slightly in error (this can be as much as 10%), then you can arrive at the true ratios by multiplying each figure by 2 and rounding off; we then get 5:2:3. Clearly, the peak at 7.3 ppm, which integrates for 5, arises from the resonance of the aromatic ring protons, and the peak at 2.0 ppm, which integrates for 3, is caused by the methyl protons. The two-proton resonance at 5.1 ppm arises from the benzyl protons. Notice then that the integrals give the simplest ratios, but not necessarily the true ratios, of the number of protons in each type.

In addition to the rising integral line, modern instruments usually give digitized numerical values for the integrals. Like the heights of the integral lines, these digitized integral values are not absolute but relative, and they should be treated as explained in the preceding paragraph. These digital values are also not exact; like the integral lines, they have the potential for a small degree of error (up to 10%). Figure 26.6 is an example of an integrated spectrum of benzyl acetate determined on a 300-MHz pulsed FT-NMR instrument. The digitized values of the integrals appear under the peaks.

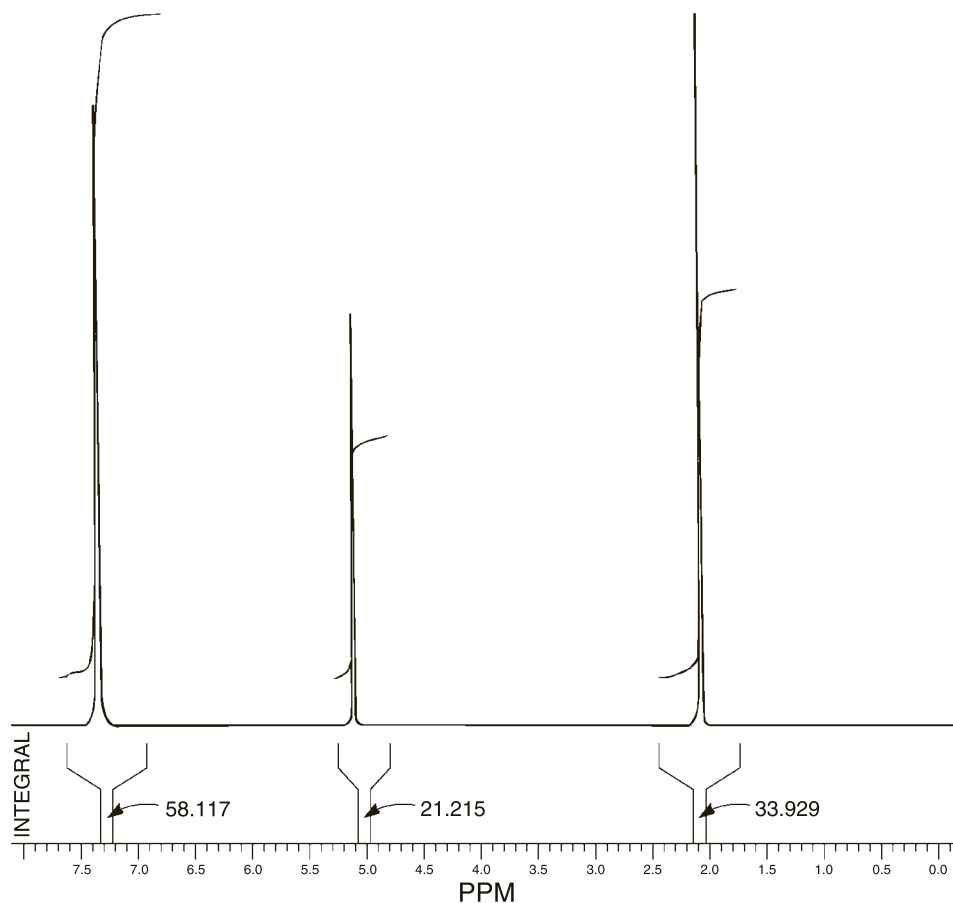


Figure 26.6

An integrated spectrum of benzyl acetate determined on a 300-MHz FT-NMR.

26.6 Chemical Environment and Chemical Shift

If the resonance frequencies of all protons in a molecule were the same, NMR would be of little use to the organic chemist. However, not only do different types of protons have different chemical shifts but they also have a value of chemical shift that characterizes the type of proton they represent. Every type of proton has only a limited range of δ values over which it gives resonance. Hence, the numerical value of the chemical shift for a proton indicates the *type of proton* originating the signal, just as the infrared frequency suggests the type of bond or functional group. Notice, for instance, that the aromatic protons of both phenylacetone (see Figure 26.3) and benzyl acetate (see Figure 26.5) have resonance near 7.3 ppm and that both methyl groups attached directly to a carbonyl group have a resonance of approximately 2.1 ppm. Aromatic protons characteristically have resonance near 7–8 ppm, and acetyl groups (the methyl protons) have their resonance near 2 ppm. These values of chemical shift are diagnostic. Notice also how the resonance of the benzyl ($-\text{CH}_2-$) protons comes at a higher value of chemical shift (5.1 ppm) in benzyl acetate than in phenylacetone (3.6 ppm). Being attached to the electronegative element, oxygen, these protons are more deshielded (see Section 26.7) than the protons in phenylacetone. A trained chemist would have readily recognized the probable presence of the oxygen by the chemical shift shown by these protons.

It is important to learn the ranges of chemical shifts over which the most common types of protons have resonance. Figure 26.7 is a correlation chart that contains the most essential and frequently encountered types of protons. Table 26.1 lists the chemical shift ranges for selected types of protons. For the beginner, it is often difficult to memorize a large body of numbers relating to chemical shifts and proton types. However, this needs to be done only approximately. It is more important to “get a feel” for the regions and the types of protons than to know a string of actual numbers. To do this, study Figure 26.7 carefully.

The values of chemical shift given in Figure 26.7 and in Table 26.1 can be easily understood in terms of two factors: local diamagnetic shielding and anisotropy. These two factors are discussed in Sections 26.7 and 26.8.

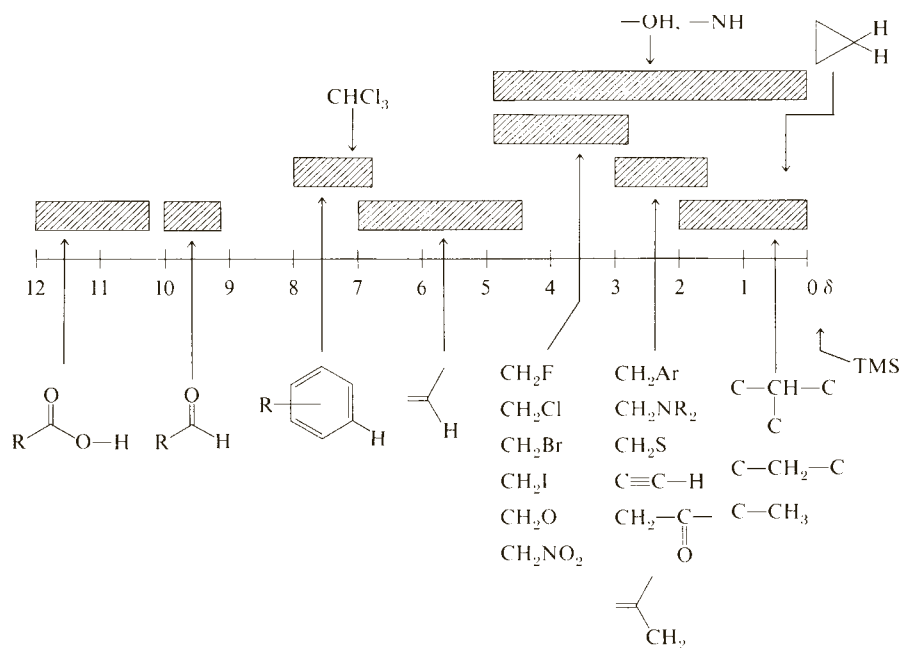
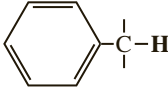
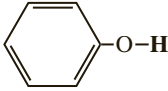
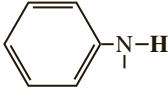
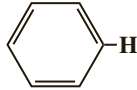


Figure 26.7

A simplified correlation chart for proton chemical shift values.

Table 26.1 Approximate Chemical Shift Ranges (ppm) for Selected Types of Protons

$R-CH_3$		0.7–1.3	$R-\overset{ }{N}-\overset{ }{C}-H$	2.2–2.9
$R-CH_2-R$		1.2–1.4	$R-S-\overset{ }{C}-H$	2.0–3.0
R_3CH		1.4–1.7		
$R-\overset{ }{C}=\overset{ }{C}-\overset{ }{C}-H$		1.6–2.6	$I-\overset{ }{C}-H$	2.0–4.0
$R-\overset{O}{\parallel}{C}-\overset{ }{C}-H, H-\overset{O}{\parallel}{C}-\overset{ }{C}-H$		2.1–2.4	$Br-\overset{ }{C}-H$	2.7–4.1
$RO-\overset{O}{\parallel}{C}-\overset{ }{C}-H, HO-\overset{O}{\parallel}{C}-\overset{ }{C}-H$		2.1–2.5	$Cl-\overset{ }{C}-H$	3.1–4.1
$N\equiv C-\overset{ }{C}-H$		2.1–3.0	$R-\overset{O}{\parallel}{S}-\overset{O}{\parallel}{O}-\overset{ }{C}-H$	ca. 3.0
		2.3–2.7	$RO-\overset{ }{C}-H, HO-\overset{ }{C}-H$	3.2–3.8
$R-C\equiv C-H$		1.7–2.7	$R-\overset{O}{\parallel}{C}-O-\overset{ }{C}-H$	3.5–4.8
$R-S-H$	var	1.0–4.0 ^a	$O_2N-\overset{ }{C}-H$	4.1–4.3
$R-N-H$	var	0.5–4.0 ^a	$F-\overset{ }{C}-H$	4.2–4.8
$R-O-H$	var	0.5–5.0 ^a		
	var	4.0–7.0 ^a	$R-\overset{ }{C}=\overset{ }{C}-H$	4.5–6.5
	var	3.0–5.0 ^a		6.5–8.0
$R-\overset{O}{\parallel}{C}-N-H$	var	5.0–9.0 ^a	$R-\overset{O}{\parallel}{C}-H$	9.0–10.0
			$R-\overset{O}{\parallel}{C}-OH$	11.0–12.0

Note: For those hydrogens shown as $-\overset{|}{C}-H$, if that hydrogen is part of a methyl group (CH_3), the shift is generally at the low end of the range given; if the hydrogen is in a methylene group ($-CH_2-$), the shift is intermediate; and if the hydrogen is in a methine group ($-CH-$), the shift is typically at the high end of the range given.

^aThe chemical shift of these groups is variable, depending on the chemical environment in the molecule and on concentration, temperature, and solvent.

Table 26.2 Dependence of Chemical Shift of CH_3X on the Element X

Compound CH_3X	CH_3F	CH_3OH	CH_3Cl	CH_3Br	CH_3I	CH_4	$(\text{CH}_3)_4\text{Si}$
Element X	F	O	Cl	Br	I	H	Si
Electronegativity of X	4.0	3.5	3.1	2.8	2.5	2.1	1.8
Chemical shift (ppm)	4.26	3.40	3.05	2.68	2.16	0.23	0

Table 26.3 Substitution Effects

	$\underline{\text{C}}\text{HCl}_3$	$\underline{\text{C}}\text{H}_2\text{Cl}_2$	$\underline{\text{C}}\text{H}_3\text{Cl}$	$-\underline{\text{C}}\text{H}_2\text{Br}$	$-\underline{\text{C}}\text{H}_2-\text{CH}_2\text{Br}$	$-\underline{\text{C}}\text{H}_2-\text{CH}_2\text{CH}_2\text{Br}$
δ (ppm)	7.27	5.30	3.05	3.3	1.69	1.25


Note: Values apply to underlined hydrogens.

26.7 Local Diamagnetic Shielding

The trend of chemical shifts that is easiest to explain is that involving electronegative elements substituted on the same carbon to which the protons of interest are attached. The chemical shift simply increases as the electronegativity of the attached element increases. This is illustrated in Table 26.2 for several compounds of the type CH_3X .

Multiple substituents have a stronger effect than a single substituent. The influence of the substituent drops off rapidly with distance. An electronegative element has little effect on protons that are more than three carbons away from it. These effects are illustrated in Table 26.3.

Electronegative substituents attached to a carbon atom, because of their electron-withdrawing effects, reduce the valence electron density around the protons attached to that carbon. These electrons *shield* the proton from the applied magnetic field. This effect, called **local diamagnetic shielding**, occurs because the applied magnetic field induces the valence electrons to circulate. This circulation generates an induced magnetic field, which *opposes* the applied field. This is illustrated in Figure 26.8. Electronegative substituents on carbon reduce the local diamagnetic shielding in the vicinity of the attached protons because they reduce the electron density around those protons. Substituents that produce this effect are said to *deshield* the proton. The greater the electronegativity of the substituent, the more the deshielding of the protons and, hence, the greater the chemical shift of those protons.

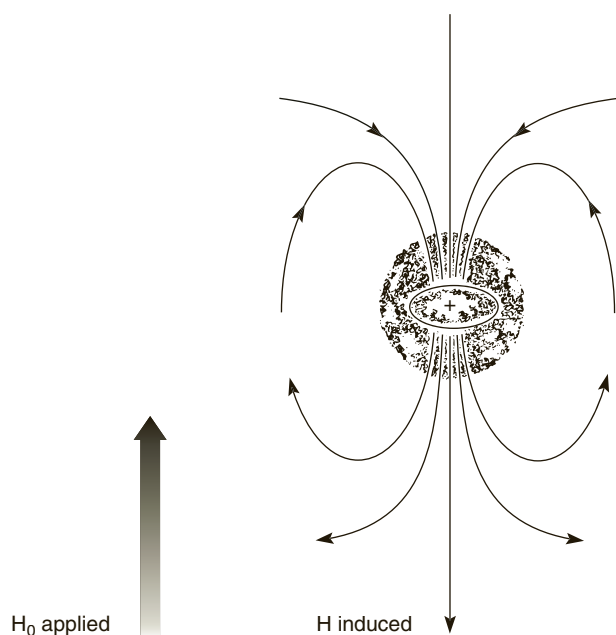**Figure 26.8**

Local diamagnetic shielding of a proton due to its valence electrons.

26.8 Anisotropy

Figure 26.7 clearly shows that several types of protons have chemical shifts not easily explained by a simple consideration of the electronegativity of the attached groups. Consider, for instance, the protons of benzene or other aromatic systems. Aryl protons generally have a chemical shift that is as large as that for the proton of chloroform. Alkenes, alkynes, and aldehydes also have protons whose resonance values are not in line with the expected magnitude of any electron-withdrawing effects. In each of these cases, the effect is due to the presence of an unsaturated

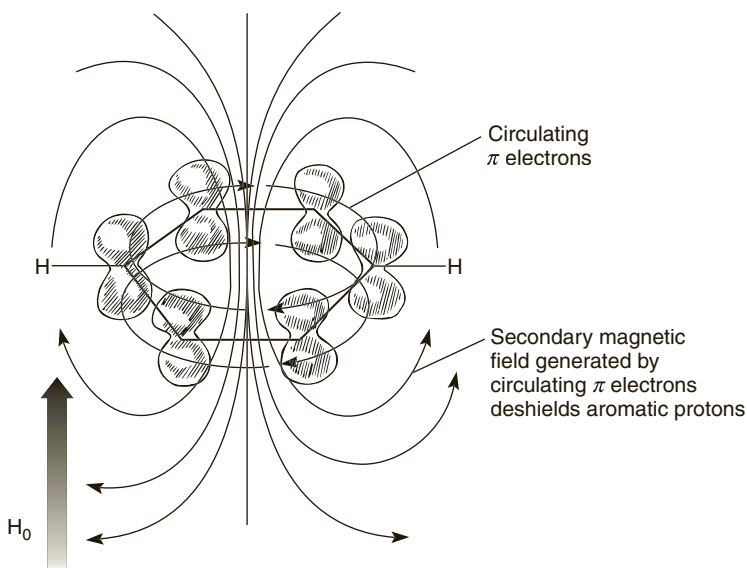


Figure 26.9
Diamagnetic anisotropy in benzene.


system (π electrons) in the vicinity of the proton in question. In benzene, for example, when the π electrons in the aromatic ring system are placed in a magnetic field, they are induced to circulate around the ring. This circulation is called a **ring current**. Moving electrons (the ring current) generate a magnetic field much like that generated in a loop of wire through which a current is induced to flow. The magnetic field covers a spatial volume large enough to influence the shielding of the benzene hydrogens. This is illustrated in Figure 26.9. The benzene hydrogens are deshielded by the **diamagnetic anisotropy** of the ring. An applied magnetic field is nonuniform (anisotropic) in the vicinity of a benzene molecule because of the labile electrons in the ring that interact with the applied field. Thus, a proton attached to a benzene ring is influenced by *three* mag-


netic fields: the strong magnetic field applied by the magnets of the NMR spectrometer and two weaker fields, one due to the usual shielding by the valence electrons around the proton and the other due to the anisotropy generated by the ring system electrons. It is this anisotropic effect that gives the benzene protons a greater chemical shift than is expected. These protons just happen to lie in a **deshielding** region of this anisotropic field. If a proton were placed in the center of the ring rather than on its periphery, the proton would be shielded because the field lines would have the opposite direction.

All groups in a molecule that have π electrons generate secondary anisotropic fields. In acetylene, the magnetic field generated by induced circulation of π electrons has a geometry such that the acetylene hydrogens are **shielded**. Hence, acetylenic hydrogens come at a higher field than expected. The shielding and deshielding regions due to the various π electron functional groups have characteristic shapes and directions; they are illustrated in Figure 26.10. Protons falling within the cones are shielded, and those falling outside the conical areas are deshielded. Because the magnitude of the anisotropic field diminishes with distance, beyond a certain distance anisotropy has essentially no effect.


26.9 Spin-Spin Splitting ($n + 1$ Rule)

We have already considered how the chemical shift and the integral (peak area) can give information about the numbers and types of hydrogens contained in a molecule. A third type of information available from the NMR spectrum is derived from spin-spin splitting. Even in simple molecules, each type of proton rarely gives a single resonance peak. For instance, in 1,1,2-trichloroethane there are two chemically distinct types of hydrogen:

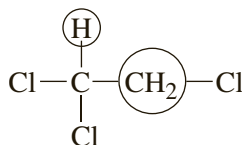

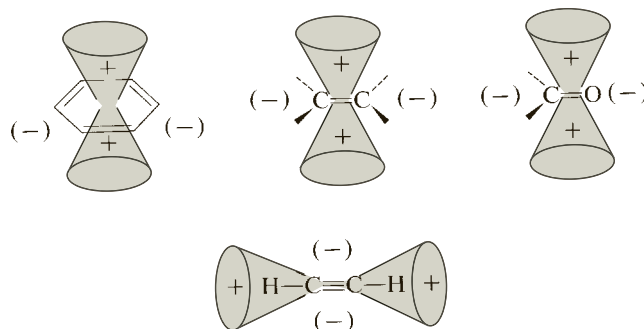**Figure 26.10**

Anisotropy caused by the presence of π electrons in some common multiple-bond systems.

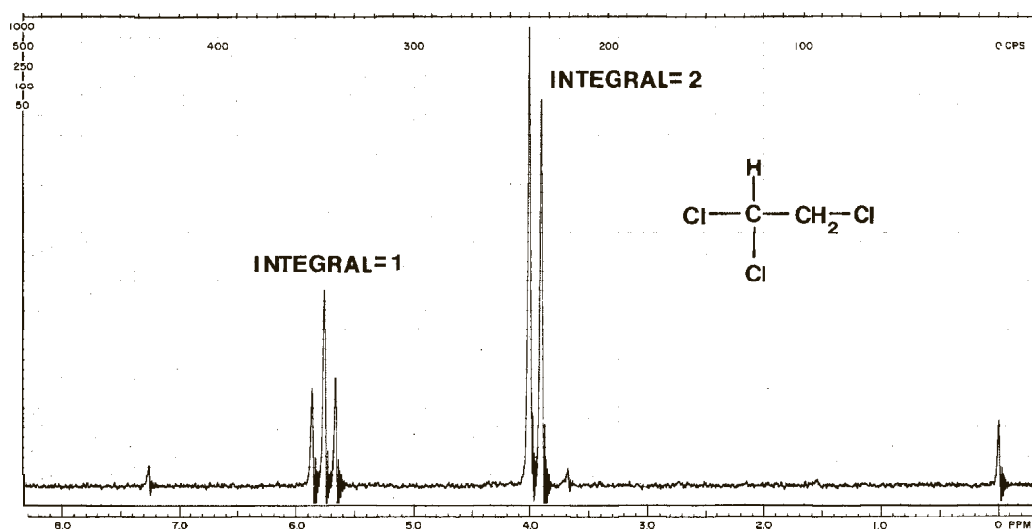**Figure 26.11**


NMR spectrum of 1,1,2-trichloroethane.

From information given thus far, you would predict *two* resonance peaks in the NMR spectrum of 1,1,2-trichloroethane with an area ratio (integral ratio) of 2:1. In fact, the NMR spectrum of this compound has *five* peaks. A group of three peaks (called a *triplet*) exists at 5.77 ppm, and a group of two peaks (called a *doublet*) is found at 3.95 ppm. The spectrum is shown in Figure 26.11. The methine (CH) resonance (5.77 ppm) is split into a triplet, and the methylene resonance (3.95 ppm) is split into a doublet. The area under the three triplet peaks is *one*, relative to an area of *two* under the two doublet peaks.

This phenomenon is called **spin-spin splitting**. Empirically, spin-spin splitting can be explained by the “ $n + 1$ rule.” Each type of proton “senses” the number of equivalent protons (n) on the carbon atom or atoms next to the one to which it is bonded, and its resonance peak is split into $n + 1$ components.

Let’s examine the case at hand, 1,1,2-trichloroethane, using the $n + 1$ rule. First, the lone methine hydrogen is situated next to a carbon bearing two methylene protons. According to the rule, it has two equivalent neighbors ($n = 2$) and is split into $n + 1 = 3$ peaks (a triplet). The methylene protons are situated next to a carbon

bearing only one methine hydrogen. According to the rule, they have one neighbor ($n = 1$) and are split into $n + 1 = 2$ peaks (a doublet).

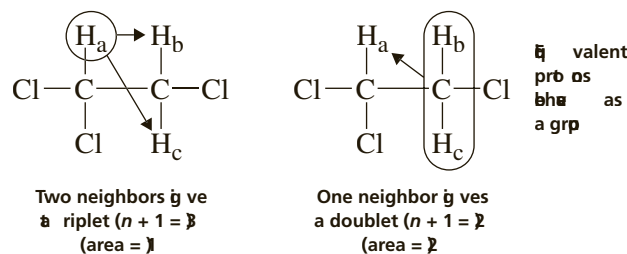

The spectrum of 1,1,2-trichloroethane can be explained easily by the interaction, or coupling, of the spins of protons on adjacent carbon atoms. The position of absorption of proton H_a is affected by the spins of protons H_b and H_c attached to the neighboring (adjacent) carbon atom. If the spins of these protons are aligned with the applied magnetic field, the small magnetic field generated by their nuclear spin properties will augment the strength of the field experienced by the first-mentioned proton H_a . The proton H_a will thus be *deshielded*. If the spins of H_b and H_c are opposed to the applied field, they will decrease the field experienced by proton H_a . It will then be *shielded*. In each of these situations, the absorption position of H_a will be altered. Among the many molecules in the solution, you will find all the various possible spin combinations for H_b and H_c ; hence, the NMR spectrum of the molecular solution will give *three* absorption peaks (a triplet) for H_a because H_b and H_c have three different possible spin combinations (Figure 26.12). By a similar analysis, it can be seen that protons H_b and H_c should appear as a doublet.

Some common splitting patterns that can be predicted by the $n + 1$ rule and that are frequently observed in a number of molecules are shown in Figure 26.13. Notice particularly the last entry, where *both* methyl groups (six protons in all) function as a unit and split the methine proton into a septet ($6 + 1 = 7$).

26.10 The Coupling Constant

The quantitative amount of spin–spin interaction between two protons can be defined by the **coupling constant**. The spacing between the component peaks in a single multiplet is called the coupling constant J . This distance is measured on the same scale as the chemical shift and is expressed in hertz (Hz).

Coupling constants for protons on adjacent carbon atoms have magnitudes of from about 6 Hz to 8 Hz (see Table 26.4). You should expect to see a coupling constant in this range for compounds where there is free rotation about a single bond. Because three bonds separate protons from each other on adjacent carbon atoms, we label these coupling constants as 3J . For example, the coupling constant for the

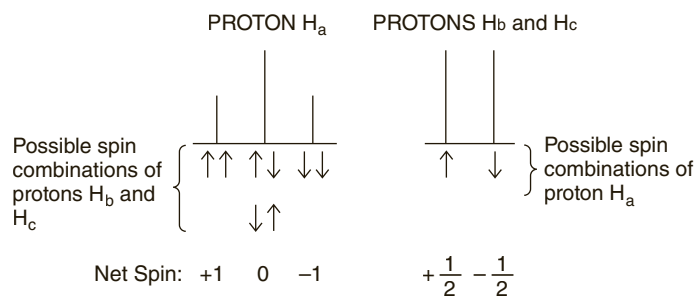


Figure 26.12
Analysis of spin–spin splitting pattern for 1,1,2-trichloroethane.

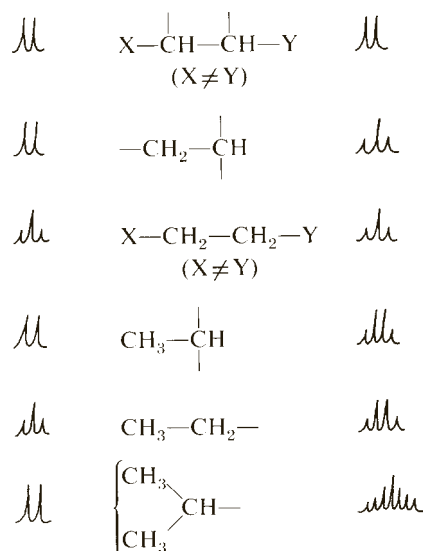



Figure 26.13
Some common splitting patterns.

compound shown in Figure 26.11 would be written as ${}^3J = 6$ Hz. The boldfaced lines in the following diagram show how the protons on adjacent carbon atoms are three bonds away from each other.


In compounds where there is a C = C double bond, free rotation is restricted. In compounds of this kind, we often find two types of 3J coupling constants: ${}^3J_{trans}$ and ${}^3J_{cis}$. These coupling constants vary in value as shown in Table 26.4, but ${}^3J_{trans}$ is

Table 26.4 Representative Coupling Constants and Approximate Values (Hz)

	3J	6–8		3J <i>ortho</i>	6–10		3J a, a	8–14
	${}^3J_{trans}$	11–18		4J <i>meta</i>	1–4		3J a, e	0–7
	${}^3J_{cis}$	6–15		5J <i>para</i>	≈ 0		3J cis	6–12
	${}^2J_{gem}$	0–5		3J <i>ortho</i>	6–10		3J trans	4–8
	3J	4–10		4J <i>meta</i>	1–4		3J cis	2–5
	3J	4–10		5J <i>para</i>	≈ 0		3J trans	1–3
	4J	0–3		5J <i>para</i>	≈ 0		3J	5–7

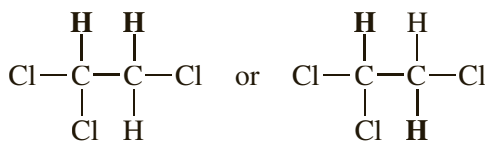

almost always larger than ${}^3J_{cis}$. The magnitudes of these 3J s often provide important structural clues. You can distinguish, for example, between a *cis* alkene and a *trans* alkene on the basis of the observed coupling constants for the two vinyl protons on disubstituted alkenes. Most of the coupling constants shown in the first column of Table 26.4 are three-bond couplings, but you will notice that there is a two-bond (2J) coupling constant listed. These protons that are bonded to a common carbon atom are often referred to as *geminal* protons and can be labeled as ${}^2J_{gem}$. Notice that the coupling constants for *geminal* protons are quite small for alkenes. The 2J couplings are observed only when the protons on a methylene group are in a different environment (see Section 26.11). The following structure shows the various types of couplings that you observe for protons on a C = C double bond in a typical alkene, vinyl acetate. The spectrum for this compound is described in detail in Section 26.11.

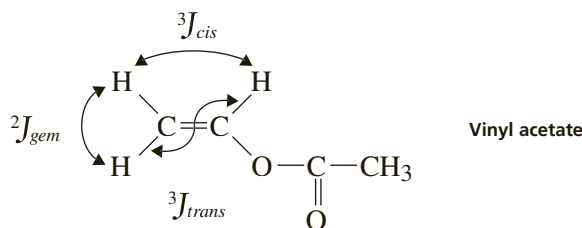


Longer-range couplings that occur over four or more bonds are observed in some alkenes and also in aromatic compounds. Thus, in Table 26.4, we see that it is possible to observe a small H—H coupling (${}^4J = 0\text{--}3$ Hz) occurring over four bonds in an alkene. In an aromatic compound, you often observe a small but measurable coupling between *meta* protons that are four bonds away from each other (${}^4J = 1\text{--}4$ Hz). Couplings over five bonds are usually quite small, with values close to 0 Hz. The long-range couplings are usually observed only in *unsaturated* compounds. The spectra of saturated compounds are often more easily interpreted because they usually have only three bond couplings. Aromatic compounds are discussed in detail in Section 26.13.

26.11 Magnetic Equivalence

In the example of spin–spin splitting in 1,1,2-trichloroethane (Figure 26.11), notice that the two protons H_b and H_c , which are attached to the same carbon atom, do not split one another. They behave as an integral group. Actually, the two protons H_b and H_c are coupled to one another; however, for reasons we cannot explain fully here, protons that are attached to the same carbon and both of which have the same chemical shift do not show spin–spin splitting. Another way of stating this is that protons coupled to the same extent to *all* other protons in a molecule do not show spin–spin splitting. Protons that have the same chemical shift and are coupled equivalently to all other protons are magnetically equivalent and do not show spin–spin splitting. Thus, in 1,1,2-trichloroethane (see Figure 26.11), protons H_b and H_c have the same value of δ and are coupled by the same value of J to proton H_a . They are magnetically equivalent, and ${}^2J_{gem} = 0$.

It is important to differentiate magnetic equivalence and chemical equivalence. Note the following two compounds.

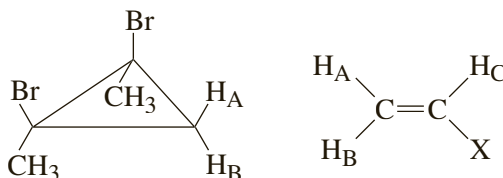

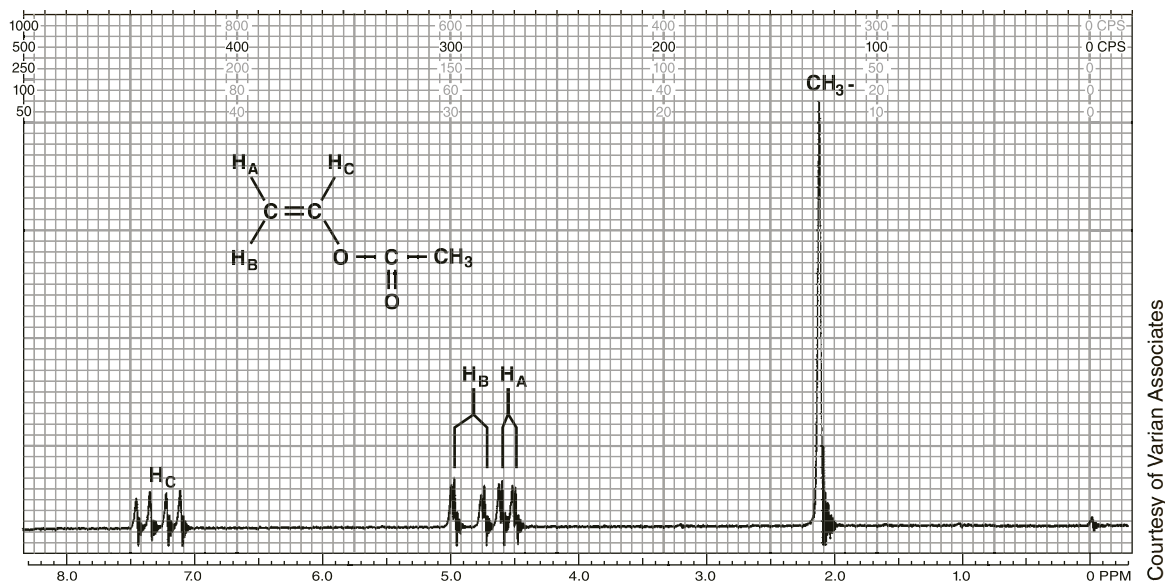


Figure 26.14
NMR spectrum of vinyl acetate.

In the cyclopropane compound, the two geminal hydrogens H_A and H_B are chemically equivalent; however, they are not magnetically equivalent. Proton H_A is on the same side of the ring as the two halogens. Proton H_B is on the same side of the ring as the two methyl groups. Protons H_A and H_B will have different chemical shifts, will couple to one another, and will show spin–spin splitting. Two doublets will be seen for H_A and H_B . For cyclopropane rings, ${}^2J_{gem}$ is usually around 5 Hz.

The general vinyl structure (alkene) shown in the previous figure and the specific example of vinyl acetate shown in Figure 26.14 are examples of cases in which the methylene protons H_A and H_B are nonequivalent. They appear at different chemical shift values and will split each other. This coupling constant, ${}^2J_{gem}$, is usually small with vinyl compounds (about 2 Hz).

The spectrum of vinyl acetate is shown in Figure 26.14. H_C appears downfield at about 7.3 ppm because of the electronegativity of the attached oxygen atom. This

proton is split by H_B into a doublet (${}^3J_{trans} = {}^3J_{BC} = 15$ Hz), and then each leg of the doublet is split by H_A into a doublet (${}^3J_{cis} = {}^3J_{AC} = 7$ Hz). Notice that the $n + 1$ rule is applied individually to each adjacent proton. The pattern that results is usually referred to as a doublet of doublets (dd). The graphic analysis shown in Figure 26.15 should help you understand the pattern obtained for proton H_C .

Now look at the pattern shown in Figure 26.14 for proton H_B at 4.85 ppm. It is also a doublet of doublets. Proton H_B is split by proton H_C into a doublet (${}^3J_{trans} = {}^3J_{BC} = 15$ Hz), and then each leg of the doublet is split by the geminal proton H_A into doublets (${}^2J_{gem} = {}^2J_{AB} = 2$ Hz).

Proton H_A , shown in Figure 26.14, appears at 4.55 ppm. This pattern is also a doublet of doublets. Proton H_A is split by proton H_C into a doublet (${}^3J_{cis} = {}^3J_{AC} = 7$ Hz), and then each leg of the

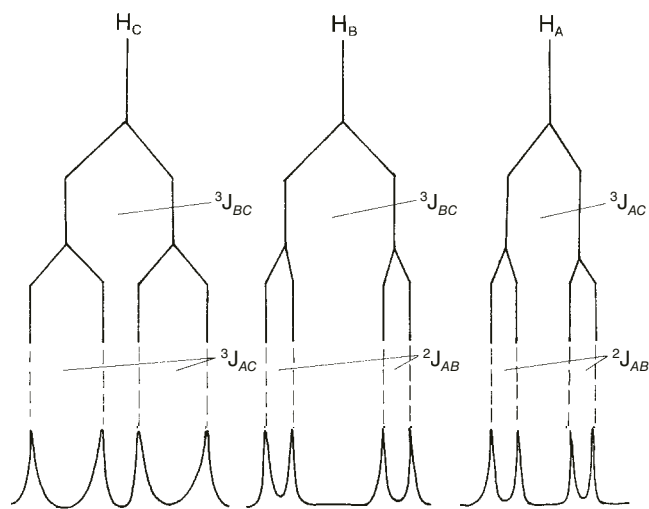


Figure 26.15
Analysis of the splittings in vinyl acetate.

doublet is split by the geminal proton H_B into doublets (${}^2J_{gem} = {}^2J_{AB} = 2$ Hz). For each proton shown in Figure 26.14, the NMR spectrum must be analyzed graphically, splitting by splitting. This complete graphic analysis is shown in Figure 26.15.

26.12 Spectra at Higher Field Strength

Occasionally, the 60-MHz spectrum of an organic compound, or a portion of it, is almost undecipherable because the chemical shifts of several groups of protons are all very similar. In these cases, all the proton resonances occur in the same area of the spectrum, and peaks often overlap so extensively that individual peaks and splittings cannot be extracted. One way to simplify such a situation is to use a spectrometer that operates at a higher frequency. Although both 60-MHz and 100-MHz instruments are still in use, it is becoming increasingly common to find instruments operating at much higher fields and with spectrometer frequencies 300, 400, or 500 MHz.

Although NMR coupling constants do not depend on the frequency or the field strength of operation of the NMR spectrometer, chemical shifts in hertz depend on these parameters. This circumstance can often be used to simplify an otherwise undecipherable spectrum. Suppose, for instance, that a compound contained three multiplets derived from groups of protons with very similar chemical shifts. At 60 MHz, these peaks might overlap, as illustrated in Figure 26.16, and simply give an unresolved envelope of absorption. It turns out that the $n + 1$ rule fails to make the proper predictions when chemical shifts are similar for the protons in a molecule. The spectral patterns that result are said to be **second order**, and what you end up seeing is an amorphous blob of unrecognizable patterns!

Figure 26.16 also shows the spectrum of the same compound at two higher frequencies (100 MHz and 300 MHz). When the spectrum is redetermined at a higher frequency, the coupling constants (J) do not change, but the chemical shifts in *hertz* (not ppm) of the proton groups (H_A , H_B , H_C) responsible for the multiplets do increase. It is important to realize, however, that the chemical shift in *ppm* is a constant, and it will not change when the frequency of the spectrometer is increased (see equation 1 in Section 26.4).

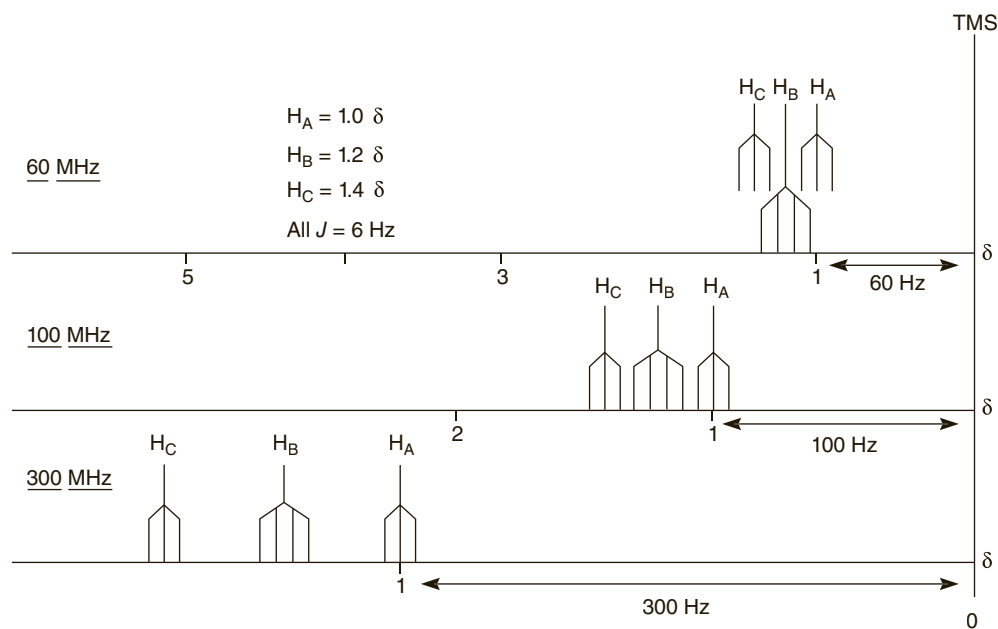


Figure 26.16

A comparison of the spectrum of a compound with overlapping multiplets at 60 MHz, with spectra of the same compound also determined at 100 MHz and 300 MHz.

Notice that at 300 MHz, the individual multiplets are cleanly separated and resolved. At high frequency, the chemical shift differences of each proton increase, resulting in more clearly recognizable patterns (that is, triplets, quartets, and so on) and less overlap of proton patterns in the spectrum. At high frequency, the chemical shift differences are large, and the $n + 1$ rule will more likely correctly predict the patterns. Thus, it is a clear advantage to use NMR spectrometers operating at high frequency (300 MHz or above) because the resulting spectra are more likely to provide nonoverlapped and well-resolved peaks. When the protons in a spectrum follow the $n + 1$ rule, the spectrum is said to be **first order**. The result is that you will obtain a spectrum with much more recognizable patterns, as shown in Figure 26.16.

26.13 Aromatic Compounds—Substituted Benzene Rings


Phenyl rings are so common in organic compounds that it is important to know a few facts about NMR absorptions in compounds that contain them. In general, the ring protons of a benzenoid system have resonance near 7.3 ppm; however, electron-withdrawing ring substituents (for example, nitro, cyano, carboxyl, or carbonyl) move the resonance of these protons downfield (larger ppm values), and electron-donating ring substituents (for example, methoxy or amino) move the resonance of these protons upfield (smaller ppm values). Table 26.5 shows these trends for a series of symmetrically *p*-disubstituted benzene compounds. The *p*-disubstituted compounds were chosen because their two planes of symmetry render all of the hydrogens equivalent. Each compound gives only one aromatic peak (a singlet) in the proton NMR spectrum. Later, you will see that some positions are affected more strongly than others in systems with substitution patterns different from this one.

In the sections that follow, we will attempt to cover some of the most important types of benzene ring substitution. In some cases, it will be necessary to examine sample spectra taken at both 60 MHz and 300 MHz. Many benzenoid rings show second-order splittings at 60 MHz, but are essentially first order at 300 MHz.

A. Monosubstituted Rings

Alkylbenzenes. In monosubstituted benzenes in which the substituent is neither a strongly electron-withdrawing nor a strongly electron-donating group, all the ring protons give rise to what appears to be a *single resonance* when the spectrum is determined at 60 MHz. This is a particularly common occurrence in alkyl-substituted benzenes. Although the protons *ortho*, *meta*, and *para* to the substituent are not chemically equivalent, they generally give rise to a single unresolved absorption peak. A possible explanation is that the chemical shift differences, which should be

Table 26.5 Proton Chemical Shifts in *p*-disubstituted Benzene Compounds

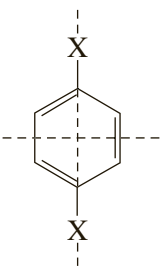Substituent X	δ (ppm)		
	—OCH ₃	6.80	} Electron-donating (shielding)
	—OH	6.60	
	—NH ₂	6.36	
	—CH ₃	7.05	
	—H	7.32	
	—COOH	8.20	} Electron-withdrawing (deshielding)
	—NO ₂	8.48	


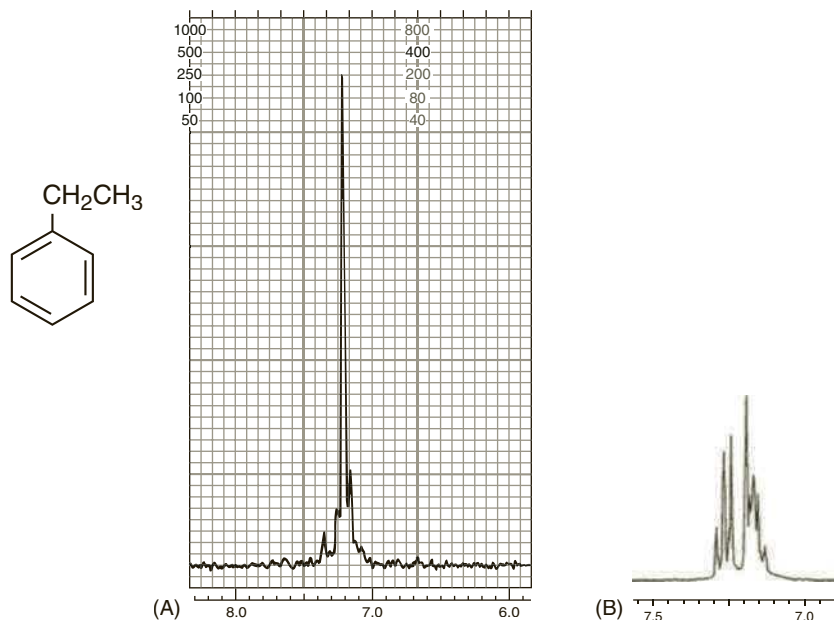


Figure 26.17
The aromatic ring portions of the ¹H NMR spectra of ethylbenzene at (A) 60 MHz and (B) 300 MHz.


small in any event, are somehow eliminated by the presence of the ring current, which tends to equalize them. All of the protons are nearly equivalent under these conditions. The NMR spectra of the aromatic portions of alkylbenzene compounds are good examples of this type of circumstance. Figure 26.17A is the 60-MHz ¹H spectrum of ethylbenzene.

The 300-MHz spectrum of ethylbenzene, shown in Figure 26.17B, presents quite a different picture. With the increased frequency shifts at 300 MHz, the nearly equivalent (at 60 MHz) protons are neatly separated into two groups. The *ortho* and *para* protons appear upfield from the *meta* protons. The splitting pattern is clearly second order.


Electron-Donating Groups. When electron-donating groups are attached to the ring, the ring protons are not equivalent, even at 60 MHz. A highly activating substituent such as methoxy clearly increases the electron density at the *ortho* and *para* positions of the ring (by resonance) and helps to give these protons greater shielding than those in the *meta* positions and, thus, a substantially different chemical shift.

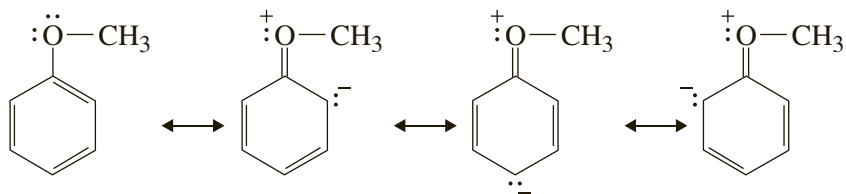

At 60 MHz, this chemical shift difference results in a complicated second-order splitting pattern for anisole (methoxybenzene), but the protons do fall clearly into two groups, the *ortho/para* protons and the *meta* protons. The 60-MHz NMR spectrum of the aromatic portion of anisole (see Figure 26.18A) has a complex multiplet for the *o,p*, protons (integrating for three protons) that is upfield from the *meta* protons

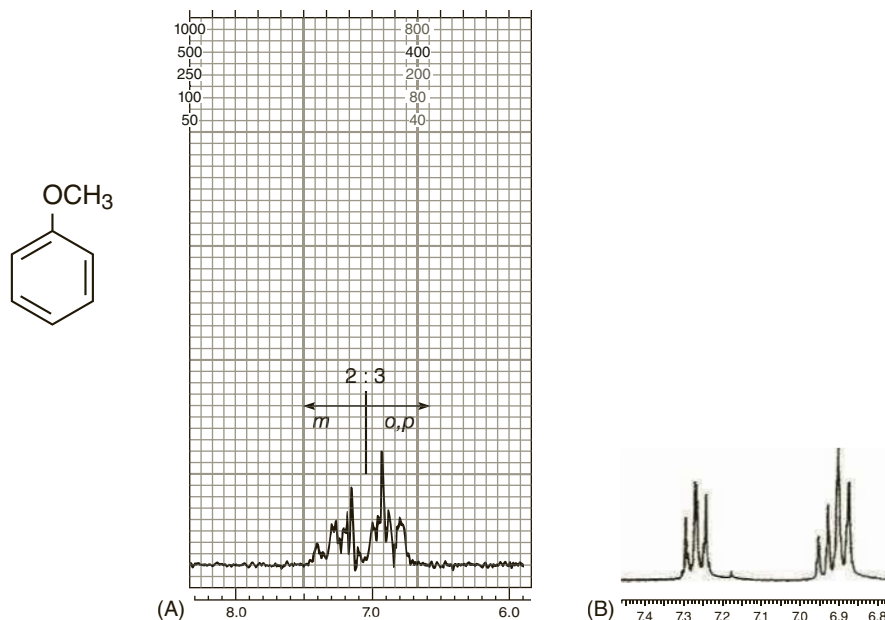**Figure 26.18**

The aromatic ring portions of the ^1H NMR spectra of anisole at (A) 60 MHz and (B) 300 MHz.

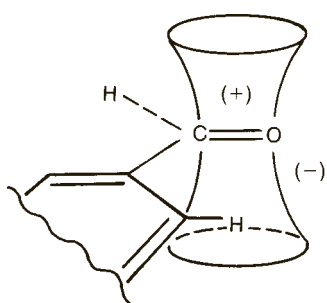**Figure 26.19**

Anisotropic deshielding of the *ortho* protons of benzaldehyde.

(integrating for two protons), with a clear distinction (gap) between the two types. Aniline (aminobenzene) provides a similar spectrum, also with a 3:2 split, owing to the electron-releasing effect of the amino group.

The 300-MHz spectrum of anisole (see Figure 26.18B) shows the same separation between the *ortho/para* hydrogens (upfield) and the *meta* hydrogens (downfield). However, because the actual shift in Hertz between the two types of hydrogens is greater, there is less second-order interaction and the lines in the pattern are sharper at 300 MHz. In fact, it might be tempting to try to interpret the observed pattern as if it were first order, a triplet at 7.25 ppm (*meta*, 2 H) and an overlapping triplet (*para*, 1 H) with a doublet (*ortho*, 2 H) at about 6.9 ppm.

Anisotropy—Electron-Withdrawing Groups. A carbonyl or a nitro group would be expected to show (aside from anisotropy effects) a reverse effect, because these groups are electron withdrawing. It would be expected that the group would act to decrease the electron density around the *ortho* and *para* positions, thus deshielding the *ortho* and *para* hydrogens and providing a pattern exactly the reverse of the one shown for anisole (3:2 ratio, downfield:upfield). Convince yourself of this by drawing resonance structures. Nevertheless, the actual NMR spectra of nitrobenzene and benzaldehyde do not have the appearances that would be predicted on the basis of resonance structures. Instead, the *ortho* protons are much more deshielded than the *meta* and *para* protons, due to the magnetic anisotropy of the π bonds in these groups.

Anisotropy is observed when a substituent group bonds a carbonyl group directly to the benzene ring (see Figure 26.19). Once again, the ring protons fall into two groups, with the *ortho* protons downfield from the *meta/para* protons. Benzaldehyde (see Figure 26.20) and acetophenone both show this effect in their NMR spectra. A similar effect is sometimes observed when a carbon–carbon double bond is attached to the ring. The 300-MHz spectrum of benzaldehyde (see Figure 26.20B) is a nearly first-order spectrum and shows a doublet (H_C , 2 H), a triplet (H_B , 1 H), and a triplet (H_A , 2 H). It can be analyzed by the $n + 1$ rule.

B. *para*-Disubstituted Rings

Of the possible substitution patterns of a benzene ring, some are easily recognized. One of these is the *para*-disubstituted benzene ring. Examine anethole (see Figure 26.21) as a first example.

On one side of the anethole ring shown in Figure 26.21, proton H_a is coupled to H_b , $^3J = 8$ Hz, resulting in a doublet at about 6.80 ppm in the spectrum. Proton H_a appears upfield (smaller ppm value) relative to H_b because of shielding by the electron-releasing effect of the methoxy group. Likewise, H_b is coupled to H_a , $^3J = 8$ Hz, producing another doublet at 7.25 ppm for this proton. Because of the plane of symmetry, both halves of the ring are equivalent. Thus, H_a and H_b on the other side of the ring also appear at 6.80 ppm and 7.25 ppm, respectively. Each doublet, therefore, integrates for two protons each. A *para*-disubstituted ring, with two different

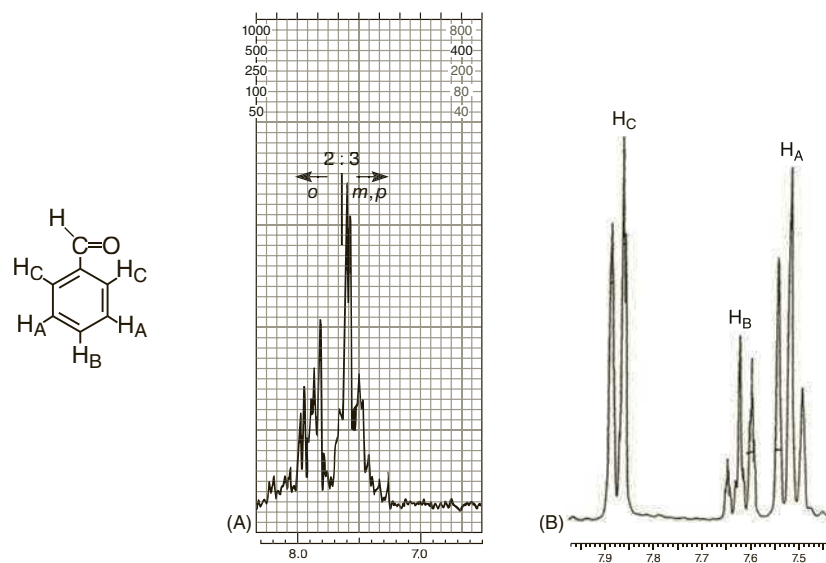


Figure 26.20

The aromatic ring portions of the ^1H NMR spectra of benzaldehyde at (A) 60 MHz and (B) 300 MHz.

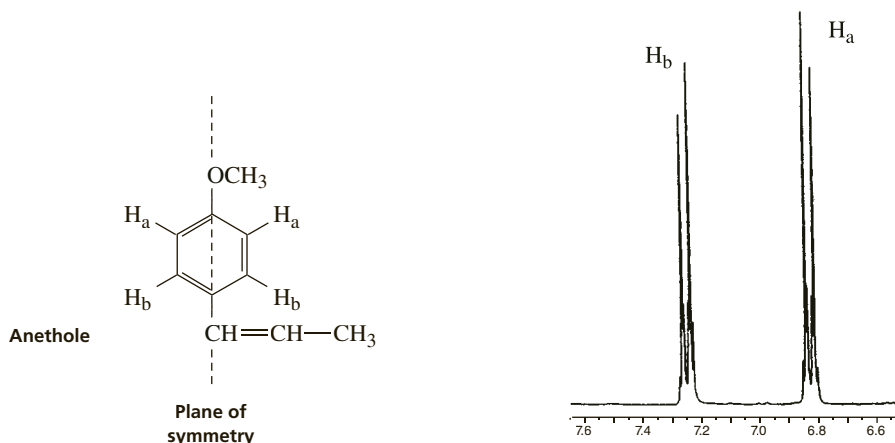
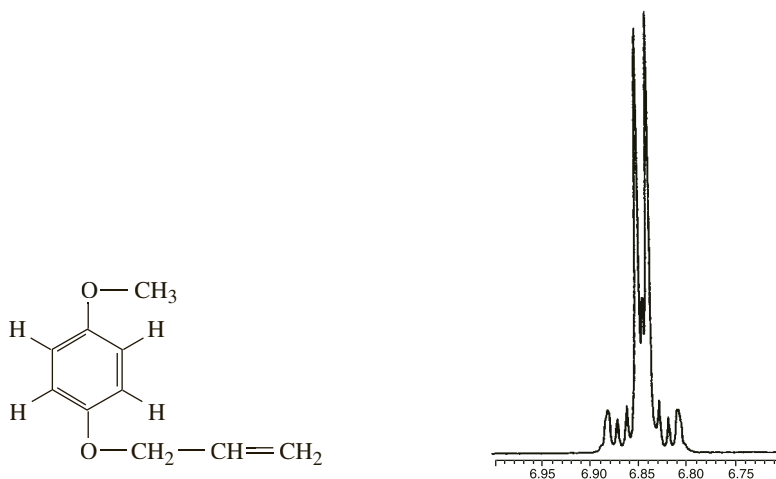


Figure 26.21

The aromatic ring protons of the 300-MHz ^1H NMR spectrum of anethole showing a *para*-disubstituted pattern.

**Figure 26.22**

The aromatic ring protons of the 300-MHz ^1H NMR spectrum of 4-allyloxyanisole.


substituents attached, is easily recognized by the appearance of two doublets, each integrating for two protons each.

As the chemical shifts of H_a and H_b approach each other in value, the *para*-disubstituted pattern becomes similar to that of 4-allyloxyanisole (see Figure 26.22). The inner peaks move closer together, and the outer ones become smaller or even disappear. Ultimately, when H_a and H_b approach each other closely enough in chemical shift, the outer peaks disappear, and the two inner peaks merge into a *singlet*; 1,4-dimethylbenzene (*para*-xylene), for instance, gives a singlet at 7.05 ppm. Hence, a single aromatic resonance integrating for four protons could easily represent a *para*-disubstituted ring, but the substituents would obviously be either identical or very similar.

C. Other Substitution

Figure 26.23 shows the 300-MHz ^1H spectra of the aromatic ring portions of 2-, 3-, and 4-nitroaniline (the *ortho*, *meta*, and *para* isomers). The characteristic pattern of a *para*-disubstituted ring, with its pair of doublets, makes it easy to recognize 4-nitroaniline. The splitting patterns for 2- and 3-nitroaniline are first order, and they can be analyzed by the $n + 1$ rule. As an exercise, see if you can analyze these patterns, assigning the multiplets to specific protons on the ring. Use the indicated multiplicities (s, d, t) and expected chemical shifts to help your assignments. Remember that the amino group releases electrons by resonance, and the nitro group shows a significant anisotropy toward *ortho* protons. You may ignore any *meta* and *para* couplings, remembering that these long-range couplings will be too small in magnitude to be observed on the scale on which these figures are presented. If the spectra were expanded, you would be able to observe 4J couplings.

The spectrum shown in Figure 26.24 is of 2-nitrophenol. It is helpful to look also at the coupling constants for the benzene ring found in Table 26.4. Because the spectrum is expanded, it is now possible to see 3J couplings (about 8 Hz), as well as 4J couplings (about 1.5 Hz). 5J couplings are not observed ($^5J = 0$). Each of the protons on this compound is assigned on the spectrum. Proton H_d appears downfield at 8.11 ppm as a doublet of doublets ($^3J_{ad} = 8$ Hz and $^4J_{cd} = 1.5$ Hz); H_c appears at 7.6 ppm as a triplet of doublets ($^3J_{ac} = ^3J_{bc} = 8$ Hz and $^4J_{cd} = 1.5$ Hz); H_b appears at 7.17 ppm as a doublet of doublets ($^3J_{bc} = 8$ Hz and $^4J_{ab} = 1.5$ Hz); and H_a

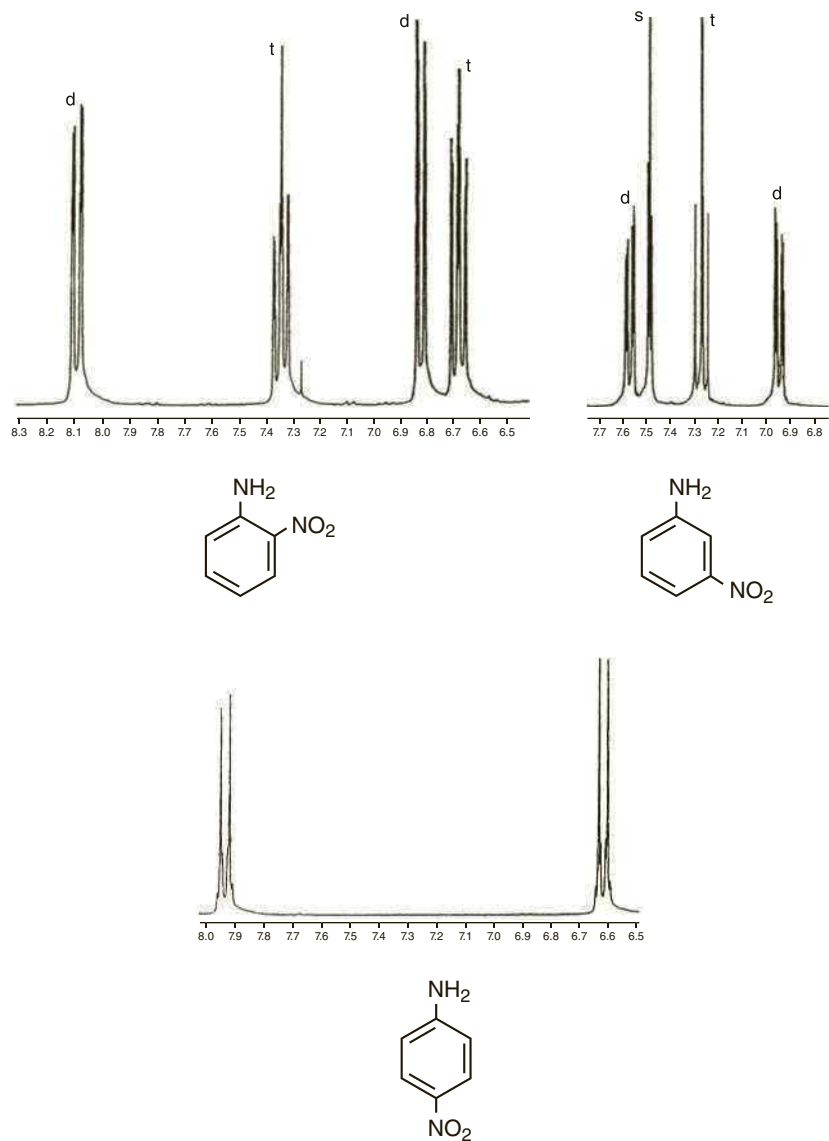
**Figure 26.23**

The 300-MHz ^1H NMR spectra of the aromatic ring portions of 2-, 3-, and 4-nitroaniline (s, singlet; d, doublet; t, triplet). The NH_2 group is not shown.

appears at 7.0 ppm as a triplet of doublets ($^3J_{ac} = ^3J_{ad} = 8$ Hz and $^4J_{ab} = 1.5$ Hz). H_d appears the furthest downfield because of the anisotropy of the nitro group. H_a and H_b are relatively shielded because of the resonance-releasing effect of the hydroxyl group, which shields these two protons. H_c is assigned by a process of elimination in the absence of these two effects.

26.14 Protons Attached to Atoms Other Than Carbon

Protons attached to atoms other than carbon often have a widely variable range of absorptions. Several of these groups are tabulated in Table 26.6. In addition, under the usual conditions of determining an NMR spectrum, protons on heteroelements normally do not couple with protons on adjacent carbon atoms to give spin–spin splitting. The primary reason is that such protons often exchange rapidly with those of the solvent medium. The absorption position is variable because these groups also undergo various degrees of hydrogen bonding in solutions of different

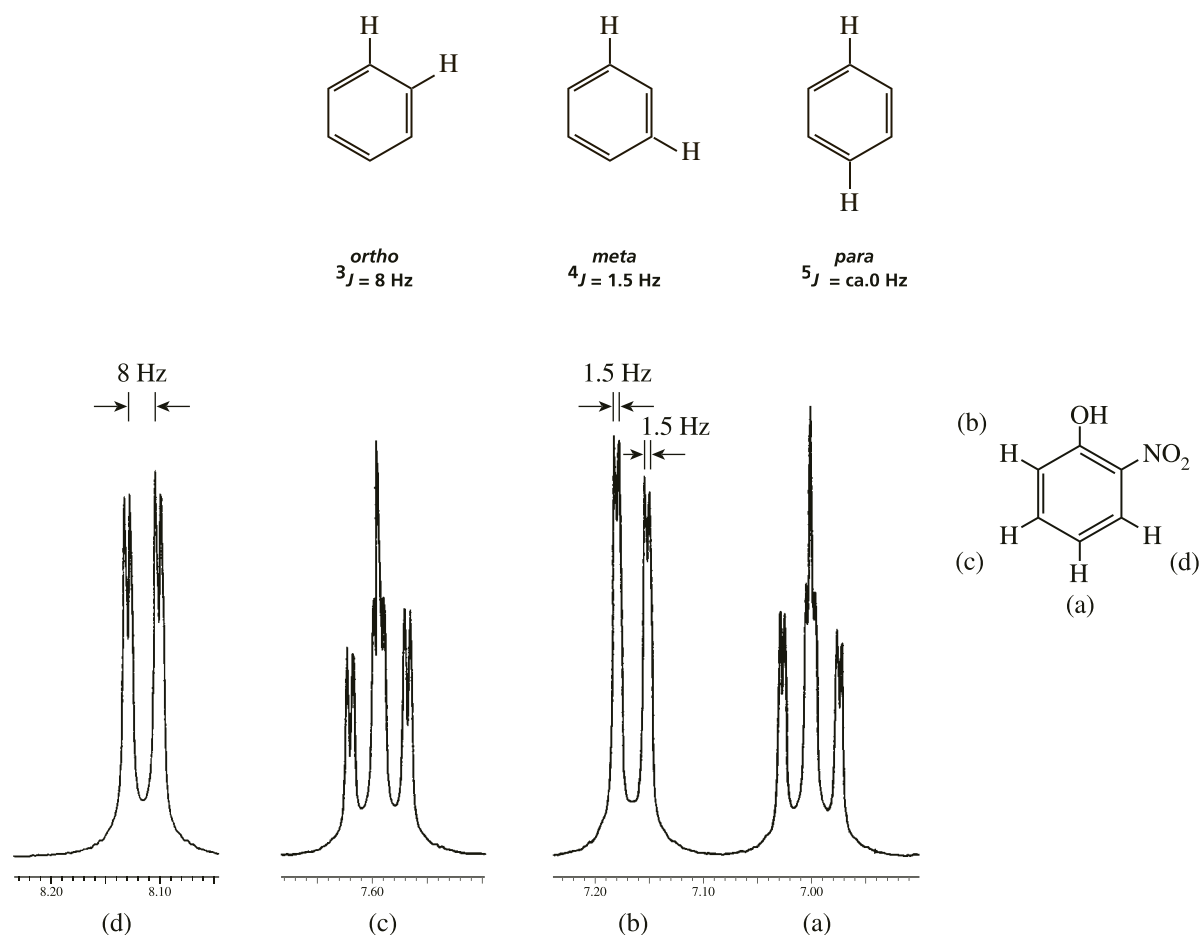**Figure 26.24**

Expansions of the aromatic ring proton multiplets from the 300-MHz ^1H spectrum of 2-nitrophenol. The accompanying hydroxyl absorption (OH) is not shown. Coupling constants are indicated on some of the peaks of the spectrum to give an idea of scale.

Table 26.6 Typical Ranges for Groups With Variable Chemical Shift


Acids	RCOOH	10.5–12.0 ppm
Phenols	ArOH	4.0–7.0
Alcohols	ROH	0.5–5.0
Amines	RNH ₂	0.5–5.0
Amides	RCONH ₂	5.0–8.0
Enols	CH=CH—OH	≥15

concentrations. The amount of hydrogen bonding that occurs with a proton radically affects the valence electron density around that proton and produces correspondingly large changes in the chemical shift. The absorption peaks for protons that have hydrogen bonding or are undergoing exchange are frequently broad relative to other singlets and can often be recognized on that basis. For a different reason, called **quadrupole broadening**, protons attached to nitrogen atoms often show an extremely broad resonance peak, often almost indistinguishable from the baseline.


26.15 Chemical Shift Reagents

Researchers have known for some time that interactions between molecules and solvents, such as those due to hydrogen bonding, can cause large changes in the resonance positions of certain types of protons (for example, hydroxyl and amino). They have also known that the resonance positions of some groups of protons can be greatly affected by changing from the usual NMR solvents such as CCl_4 and CDCl_3 to solvents such as benzene, which impose local anisotropic effects on surrounding molecules. In many cases, it is possible to resolve partially overlapping multiplets by such a solvent change. The use of **chemical shift reagents** for this purpose dates from about 1969. Most of these chemical shift reagents are organic complexes of paramagnetic rare earth metals from the lanthanide series of elements. When these metal complexes are added to the compound whose spectrum is being determined, profound shifts in the resonance positions of the various groups of protons are observed. The direction of the shift (upfield or downfield) depends primarily on which metal is being used. Complexes of europium, erbium, thulium, and ytterbium shift resonances to lower field; complexes of cerium, praseodymium, neodymium, samarium, terbium, and holmium generally shift resonances to higher field. The advantage of using such reagents is that shifts similar to those observed at higher field can be induced without the purchase of an expensive higher-field instrument.

Of the lanthanides, europium is probably the most commonly used metal. Two of its widely used complexes are *tris*-(dipivalomethanato)europium and *tris*-(6,6,7,7,8,8,8-heptafluoro-2,2-dimethyl-3,5-octanedionato)europium. These are frequently abbreviated $\text{Eu}(\text{dpm})_3$ and $\text{Eu}(\text{fod})_3$, respectively.

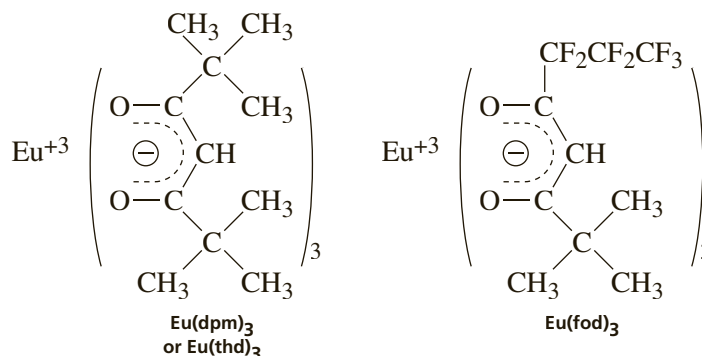


These lanthanide complexes produce spectral simplifications in the NMR spectrum of any compound that has a relatively basic pair of electrons (unshared pair) that can coordinate with Eu^{3+} . Typically, aldehydes, ketones, alcohols, thiols, ethers, and amines will all interact:

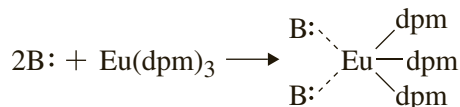


The amount of shift that a given group of protons will experience depends (1) on the distance separating the metal (Eu^{3+}) and that group of protons, and (2) on the concentration of the shift reagent in the solution. Because of the latter dependence, it is necessary when reporting a lanthanide-shifted spectrum to report the number of mole equivalents of shift reagent used or its molar concentration.

The distance factor is illustrated in the spectra of hexanol, which are given in Figures 26.25 and 26.26. In the absence of shift reagent, the normal spectrum is obtained (see Figure 26.25). Only the triplet of the terminal methyl group and the triplet of the methylene group next to the hydroxyl are resolved in the spectrum. The other protons (aside from OH) are found together in a broad unresolved group.

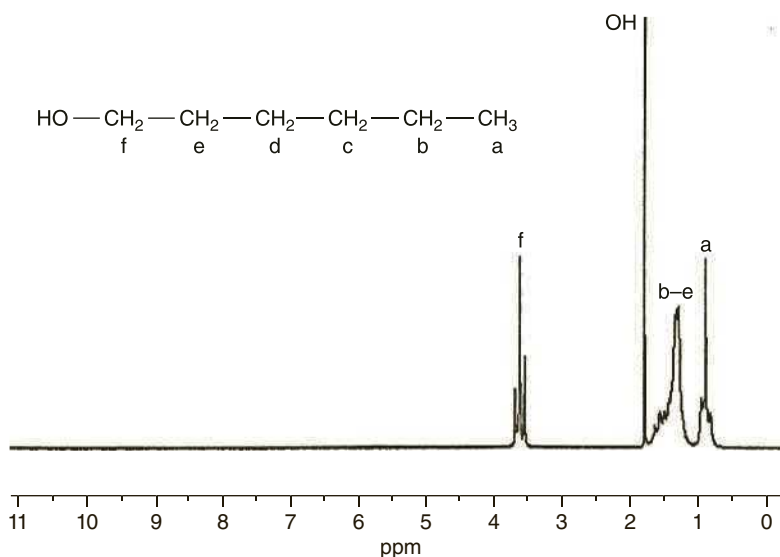**Figure 26.25**

90-MHz ^1H NMR spectrum of hexanol determined without $\text{Eu}(\text{dpm})_3$ © National Institute of Advanced Industrial Science and Technology.

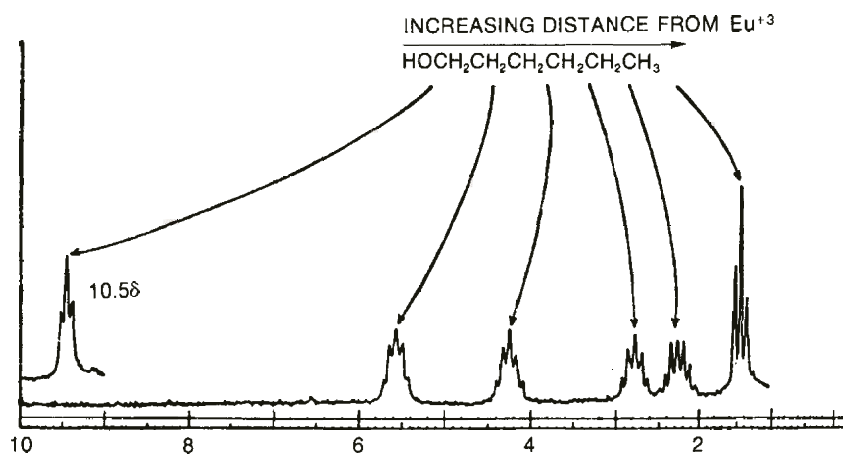**Figure 26.26**

The 100-MHz ^1H NMR spectrum of hexanol with 0.29 mole equivalents of $\text{Eu}(\text{dpm})_3$ added. From J. K. M. Sanders and D. H. Williams, *Chemical Communications*, (1970): 422. Reproduced with permission of the Royal Society of Chemistry.

With shift reagent added (see Figure 26.26), each of the methylene groups is clearly separated and resolved into the proper multiplet structure. The spectrum is first-order and simplified; all the splittings are explained by the $n + 1$ rule.

One final consequence of using a shift reagent should be noted. Notice in Figure 26.26 that the multiplets are not as nicely resolved into sharp peaks as you might expect. This is due to the fact that shift reagents cause a small amount of peak broadening. At high-shift reagent concentrations, this problem becomes serious, but at most useful concentrations the amount of broadening experienced is tolerable.

Section 26.16 Diastereopic protons

Coupling constants were discussed in Section 26.9 and 26.10 with Table 26.4 showing typical values for coupling constants in organic compounds. A simple organic compound is shown in Figure 26.27, where it is easy to predict the splitting patterns using the $n + 1$ rule. The protons H_c and H_d are equivalent, and as a group they are split by the proton H_e into a doublet, integrating for 2 H. Notice that the protons H_c and H_d are in the same electronic environment, adjacent to the two chlorine atoms, and will not split each other. *Only non-equivalent protons will split each other!* Proton H_e is split by two neighboring protons H_c and H_d into a triplet, integrating for 1 H. All of these couplings involve protons that are 3 bonds away from each other, and are referred to as 3J couplings. Similarly, the methyl group H_a appears as a triplet, integrating for 3 H, while its neighboring methylene group H_b appears as a quartet, integrating for 2 H. Many organic compounds show these simple splitting patterns. These common organic compounds show couplings between protons that are 3 bonds away from each other (3J) making the application of the $n + 1$ rule possible (see Section 26.9).

Figure 26.28 shows the structure of the product from the chiral reduction of ethyl acetoacetate (Experiment 29) the 300 MHz NMR spectrum is shown in Experiment 29. This high field NMR spectrum makes it possible to do a detailed analysis of the compound.³ The protons H_c and H_d are in different environments; H_c is opposite the OH group while H_d is further away from the OH group. The carbon atom labeled with an asterisk has a stereocenter. Often the presence of a stereocenter

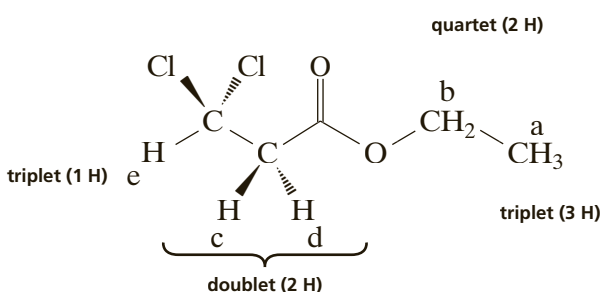


Figure 26.27

Analysis of the NMR pattern for ethyl 3,3-dichloropropanoate. Proton H_c and H_d are equivalent and do not split each other. *These methylene protons are not diastereotopic.*

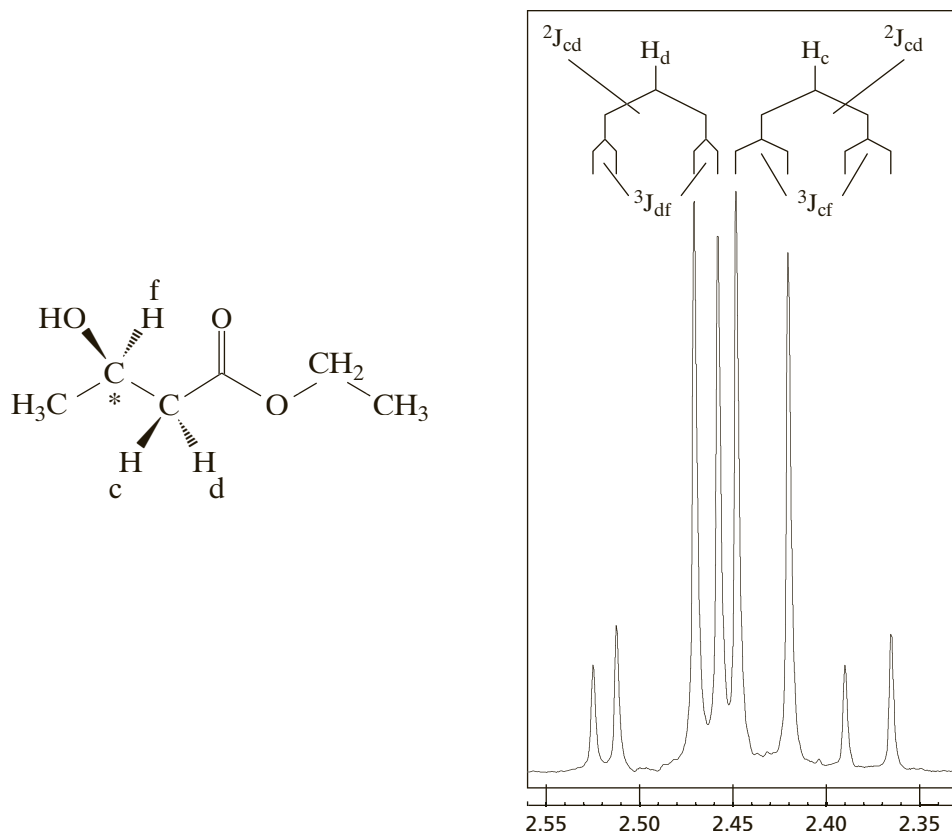
will make the protons on a methylene group non-equivalent. These non-equivalent protons are referred to as diastereotopic protons. The diastereotopic protons in compounds attached to an sp^3 hybridized carbon atom possess large coupling constants, with 2J values of 11 to 15 Hz.⁴ In fact, these values are typically larger than 3J coupling constants in compounds with sp^3 hybridization. Referring to the 300 MHz NMR spectrum of the racemic product in Experiment 29, H_c is split into a doublet by H_d and each leg of the doublet is split again into doublets by the neighboring proton, H_e . The second coupling is a three bond coupling, $^3J_{cf}$. The observed pattern is referred to as a doublet of doublets, shown in the analysis in Figure 26.28. Likewise, H_d also becomes a doublet of doublets. Both H_c and H_d appear at different places in the spectrum (See

Figure 1 and Figure 2 in Experiment 29). One may notice in Figure 26.28 that the spacing of the doublets, $^3J_{cf}$, differs from the spacing, $^3J_{df}$. The reason for these differences is that the dihedral angle for H_c to H_e is 180 degrees, while the dihedral angle for H_d to H_e is 60 degrees. The dihedral angles often influence the value of the 3J coupling constants.


Another example of a compound with diastereotopic methylene protons is shown in Figure 26.29. Notice that the hydrogen atoms on the methylene

³ Often lower field NMR spectra determined at 60 MHz are more difficult to analyze because of more closely spaced peaks.

⁴ The geminal protons (2J) attached to the terminal carbon atom of an alkene are non-equivalent and will split each other. These 2J protons attached to an sp^2 carbon atom have much smaller coupling constants, often less than 5 Hz, compared to 2J protons with sp^3 hybridized carbon atoms. (See Section 26.10 and 26.11.)

**Figure 26.28**

The structure of ethyl 3-hydroxybutanoate obtained from the reduction of ethyl 3-hydroxybutanoate (Experiment 29). This compound possesses diastereotopic methylene protons, H_c and H_d. An analysis of the diastereotopic protons is also shown in this figure.

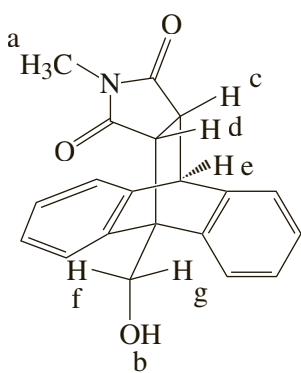**Figure 26.29**

The Diels–Alder reaction with anthracene-9-methanol from Experiment 54. This compound possesses diastereotopic methylene protons, H_f and H_g.

group are in a different environment. Proton H_f is on the same side as the nitrogen-containing ring, while H_g is on the opposite side of the nitrogen-containing ring. One observes that the molecule is not symmetrical, leading to the methylene group having diastereotopic protons. If the nitrogen-containing ring were not present, the molecule would have been symmetrical, and the methylene protons would have been equivalent and thus not diastereotopic.

This compound is produced in Experiment 54. The 500 MHz NMR spectrum of this product is shown in Figure 26.30, with insets showing expansions of the regions from 3.25 to 3.40 ppm and 4.7 to 5.2 ppm. The Hz values that appear above the peaks in the expansions may be used to calculate the coupling constants. The diastereotopic protons H_f and H_g appear at 4.96 ppm and 5.13 ppm (we don't know the assignments precisely). Assuming the correct assignment, proton H_g appears at 4.96 ppm as a doublet ($^2J_{fg} = 2490.12 - 2478.40 \text{ Hz} = \text{about } 11.7 \text{ Hz}$). However, each leg of the doublet is split again into doublets, ($^3J_{bg} = 2490.12 - 2485.36 \text{ Hz} = \text{about } 4.8 \text{ Hz}$). The 4.8 Hz coupling constant is for the 3J coupling of H_g to the proton on the OH group, H_b. The resulting pattern at 4.96 ppm is a doublet of doublets. Another doublet of doublets appears at 5.13 ppm for the non-equivalent diastereotopic proton H_f.

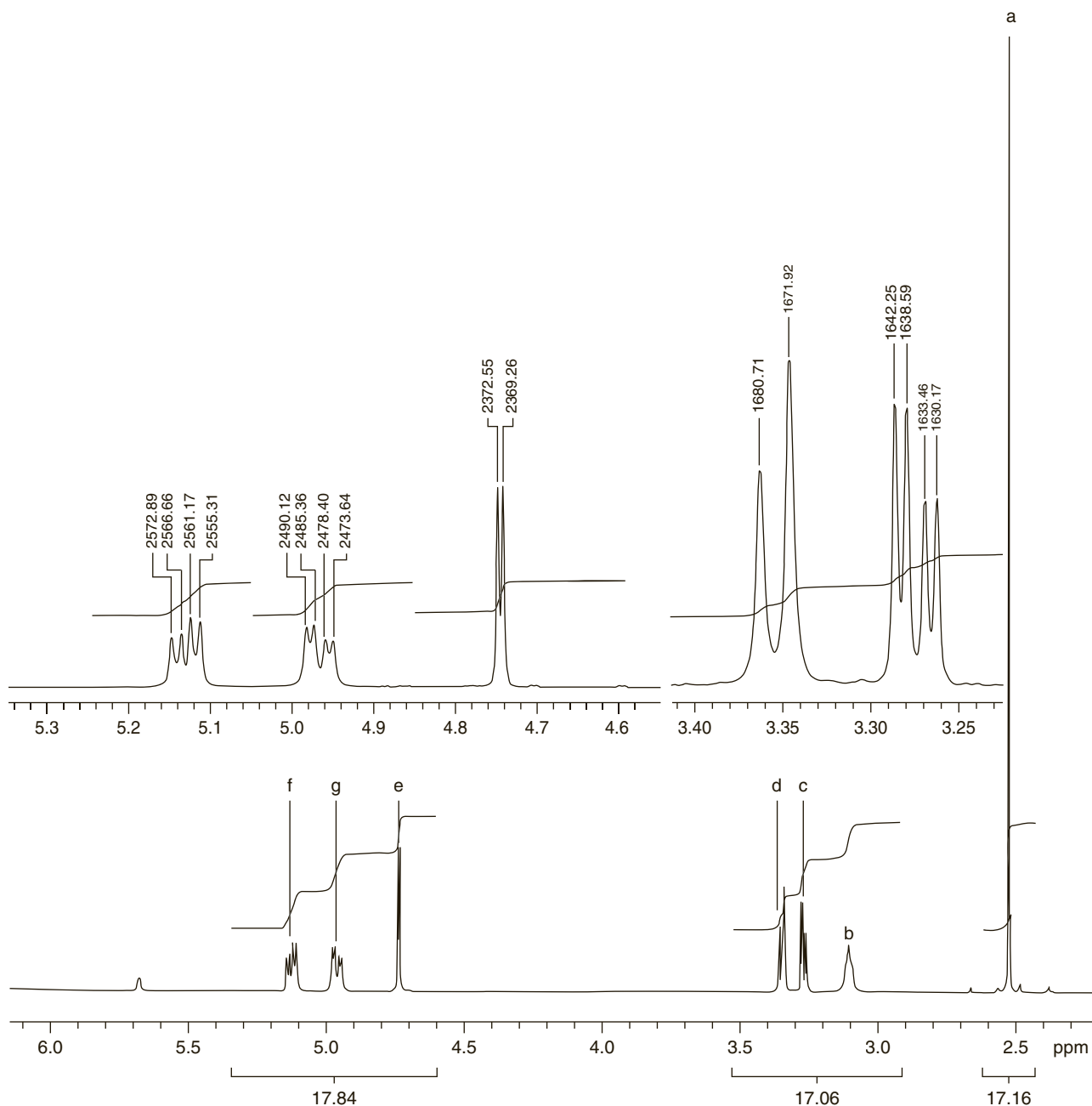


Figure 26.30

The 500 MHz NMR spectrum of the the Diels-Alder reaction with anthracene-9-methanol from Experiment 54. The inserts show expansions of the regions from 3.25 to 3.40 ppm and 4.7 to 5.2 ppm. *The eight aromatic ring protons are not shown.*

The other coupling constants may be extracted from the pattern at 3.27 ppm, a doublet of doublets, ${}^3J_{cd} = 8.8$ Hz and ${}^3J_{ce} = 3.3$ Hz. Proton H_d appears at 3.36 ppm as doublet $J_{cd} = 8.8$ Hz. Finally, H_e appears as a doublet at 4.75 ppm with ${}^3J_{ce} = 3.3$ Hz. The coupling constants may be summarized as follows:

$${}^2J_{fg} = 11.7 \text{ Hz (diastereotopic protons at 4.96 ppm and 5.13 ppm)}$$

$${}^3J_{bf} = {}^3J_{bg} = 4.8 \text{ Hz (coupling of OH to each of the diastereotopic protons)}$$

$${}^3J_{ce} = 3.3 \text{ Hz (coupling of H}_c \text{ to H}_e)$$

$${}^3J_{cd} = 8.8 \text{ Hz (coupling of H}_e \text{ to H}_d)$$

REFERENCES

Textbooks

- Friebolin, H. *Basic One- and Two-Dimensional NMR Spectroscopy*, 5th ed. VCH Publishers: New York, 2010.
- Gunther, H. *NMR Spectroscopy*, 2nd ed. John Wiley & Sons: New York, 1995.
- Jackman, L. M.; Sternhell, S. *Nuclear Magnetic Resonance Spectroscopy in Organic Chemistry*, 2nd ed. Pergamon Press: New York, 1969.
- Macomber, R. S. *A Complete Introduction to Modern NMR Spectroscopy*. John Wiley & Sons: New York, 2008.
- Mocomber, R. S. *NMR Spectroscopy: Essential Theory and Practice*. College Outline Series, Harcourt Brace Jovanovich: New York, 1998.
- Pavia, D. L.; Lampman, G. M.; Kriz, G. S.; Vyvyan, J. R. *Introduction to Spectroscopy*, 5th ed. Cengage Learning: Stamford, CT, 2015.
- Sanders, J. K.; Hunter, B. K. *Modern NMR Spectroscopy—A Guide for Chemists*, 2nd ed. Oxford University Press: Oxford, 1993.
- Silverstein, R. M.; Webster, F. X.; Kiemle, D.; Bryce, D. L. *Spectrometric Identification of Organic Compounds*, 8th ed. John Wiley & Sons: New York, 2015.

Compilations of Spectra

- Pouchert, C. J. *The Aldrich Library of NMR Spectra*, 60 MHz, 2nd ed. Aldrich Chemical Company: Milwaukee, WI, 1983.
- Pouchert, C. J.; Behnke, J. *The Aldrich Library of ^{13}C and ^1H FT-NMR Spectra*, 300 MHz. Aldrich Chemical Company: Milwaukee, WI, 1993.
- Pretsch E.; Clerc, T.; Seibl, J.; Simon, W. *Tables of Spectra Data for Structure Determination of Organic Compounds*, 2nd ed. Springer-Verlag: Berlin and New York, 1989.

Translated from the German by K. Biemann.


Web Sites

- <http://www.aist.go.jp/RIODB/SDBS/menu-e.html>
Integrated Spectral DataBase System for Organic Compounds, National Institute of Materials and Chemical Research, Tsukuba, Ibaraki 305-8565, Japan. This database includes infrared, mass spectra, and NMR data (proton and carbon-13) for a large number of compounds.
- <http://www.chem.ucla.edu/~webspectra/>
UCLA Department of Chemistry and Biochemistry in connection with Cambridge University Isotope Laboratories maintains a Web site, WebSpectra, that provides NMR and IR spectroscopy problems for students to interpret. They provide links to other sites with problems for students to solve.


PROBLEMS

1. Describe the method that you should use to determine the proton NMR spectrum of a carboxylic acid, which is insoluble in all the common organic solvents that your instructor is likely to make available.
2. To save money, a student uses chloroform instead of deuterated chloroform to run a proton NMR spectrum. Is this a good idea?
3. Look up the solubilities for the following compounds and decide whether you would select deuterated chloroform or deuterated water to dissolve the substances for NMR spectroscopy.


- a. Glycerol (1,2,3-propanetriol)
 - b. 1,4-Diethoxybenzene
 - c. Propyl pentanoate (propyl ester of pentanoic acid)
4. Assign each of the proton patterns in the spectra of 2-, 3-, and 4-nitroaniline as shown in Figure 26.23.
 5. The following two compounds are isomeric esters derived from acetic acid, each with formula $C_5H_{10}O_2$. These expanded spectra clearly show the splitting patterns: singlet, doublet, triplet, quartet, etc. Integral curves are drawn on the spectra, along with relative integration values provided just above the scale and under each set of peaks. These numbers indicate the number of protons assigned to each pattern. Remember that these integral values are approximate. You will need to round the values off to the nearest whole number. Draw the structure of each compound.

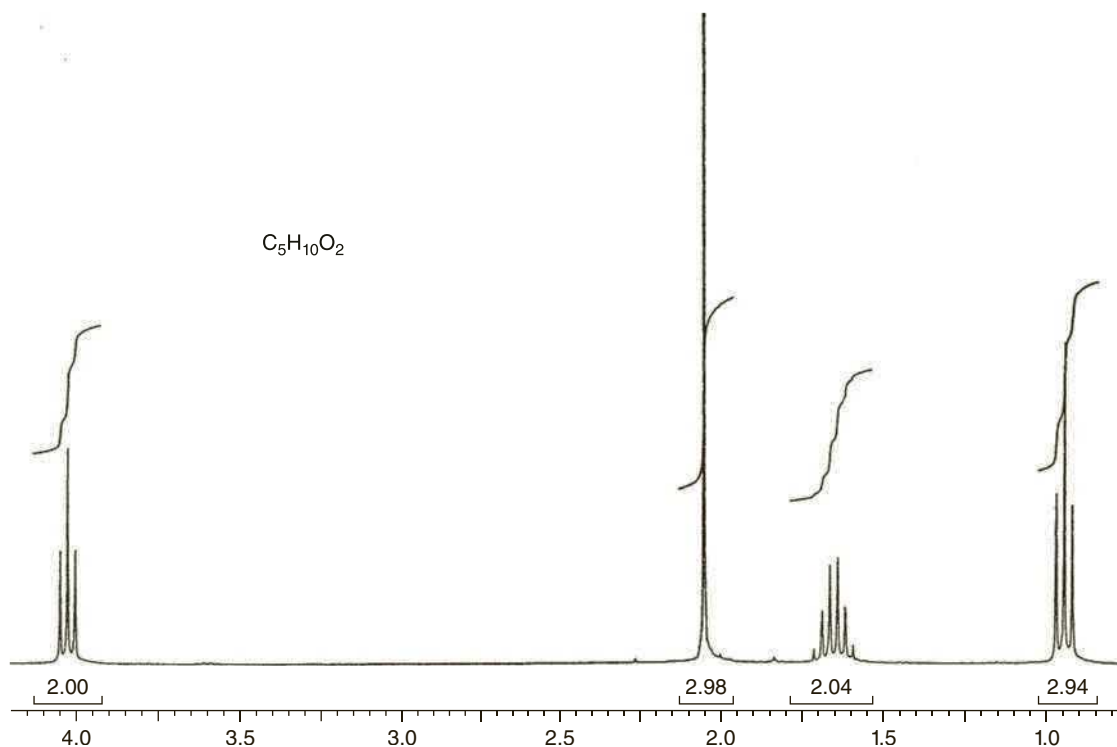


- b. The set of peaks centering on 5 ppm is expanded in both the x and y directions in order to show the pattern more clearly. This expanded pattern is shown as an inset on the full spectrum.

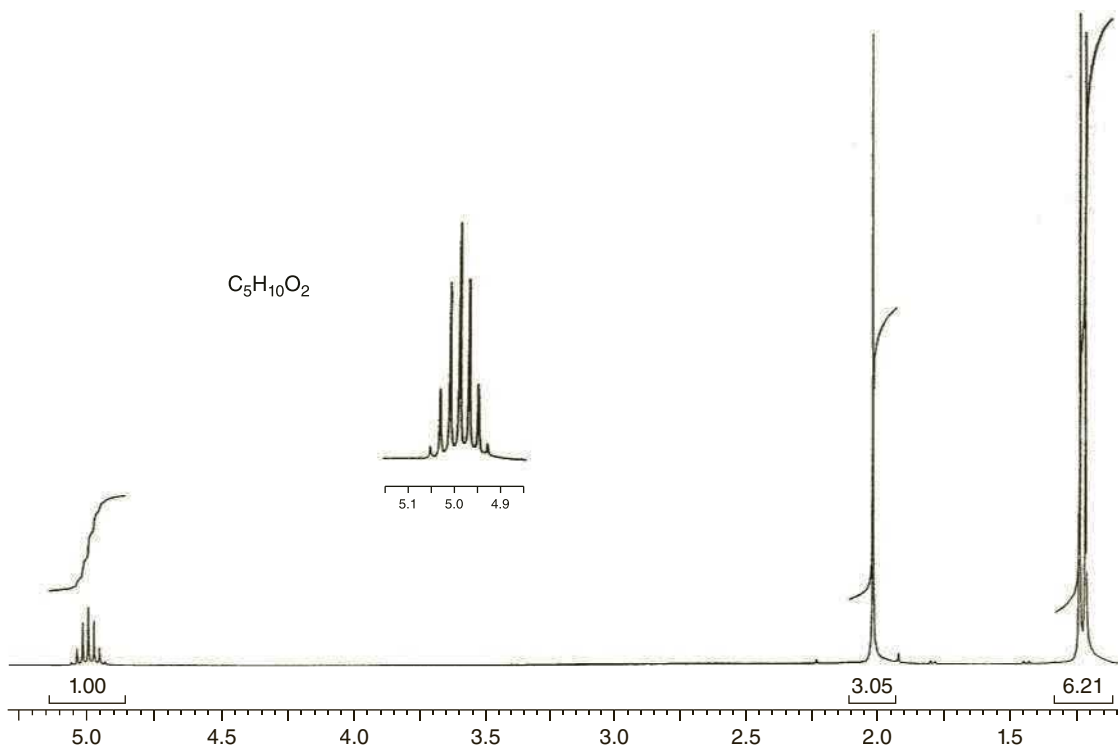

6. The compound that gives the following NMR spectrum has the formula $C_3H_6Br_2$. Draw the structure.

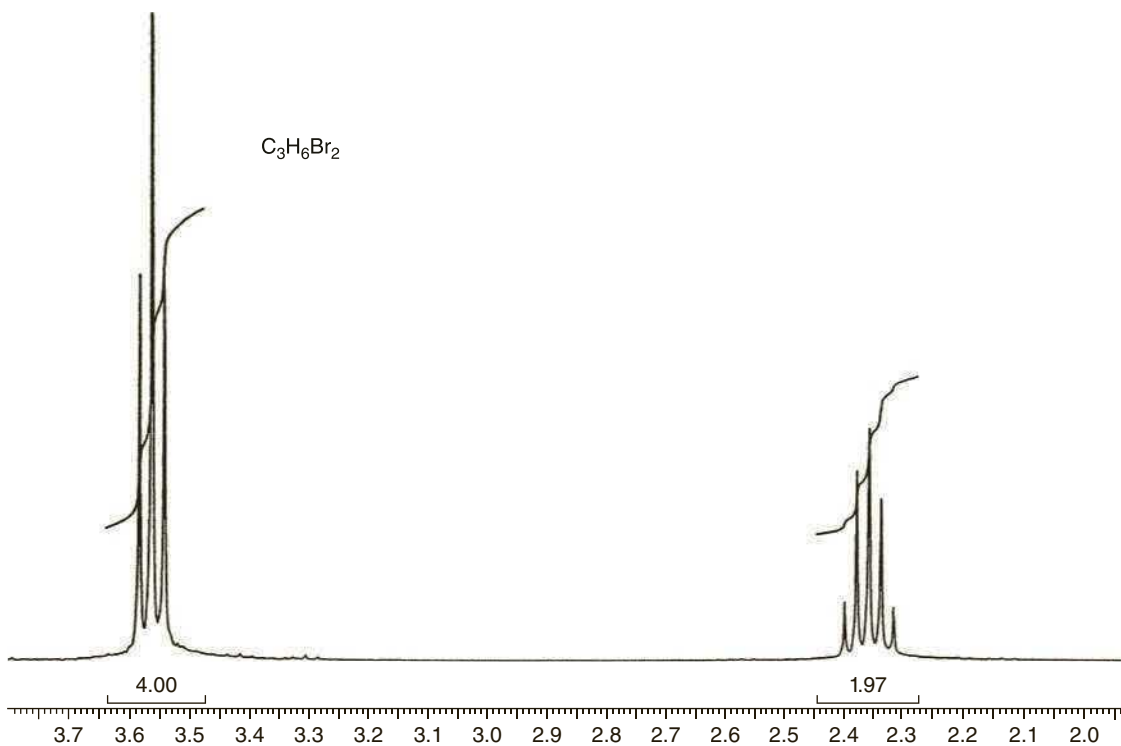




7. Draw the structure of an ether with formula $C_5H_{12}O_2$ that fits the following NMR spectrum.

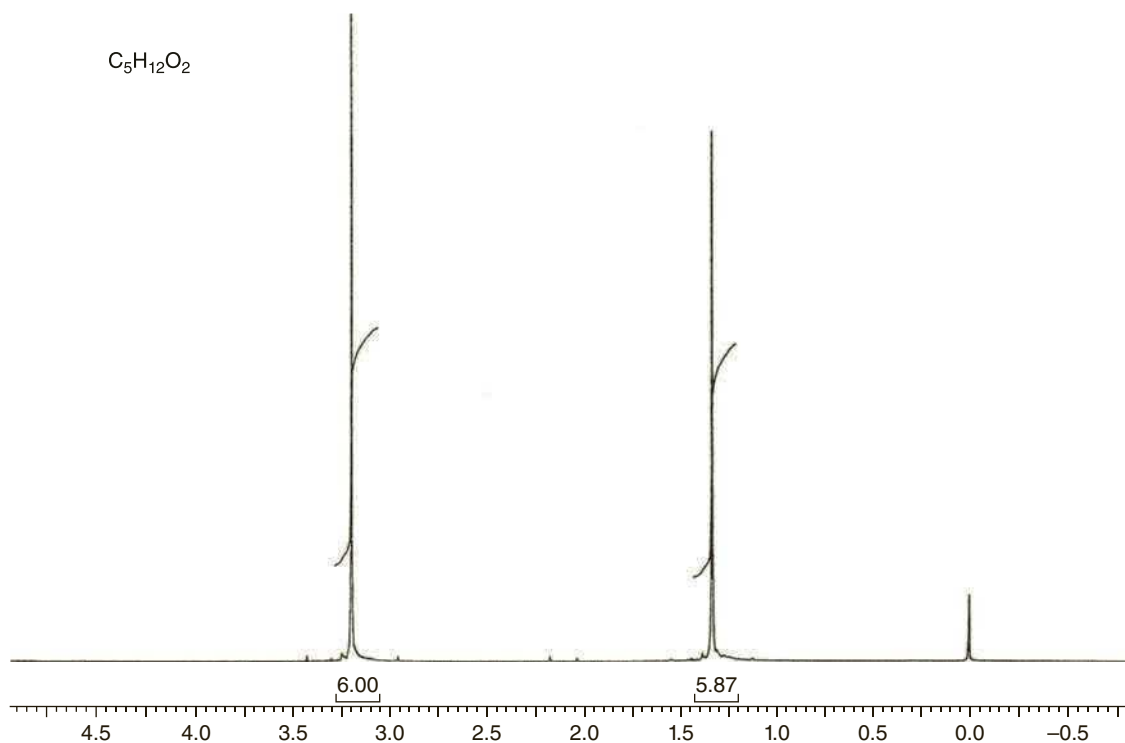

8. Following are the NMR spectra of three isomeric esters with the formula $C_7H_{14}O_2$, all derived from propanoic acid. Provide a structure for each.


- a. The set of peaks centering on about 1.9 ppm is expanded in both the x and y directions in order to show the pattern more clearly. This expanded pattern is shown as an inset on the full spectrum.

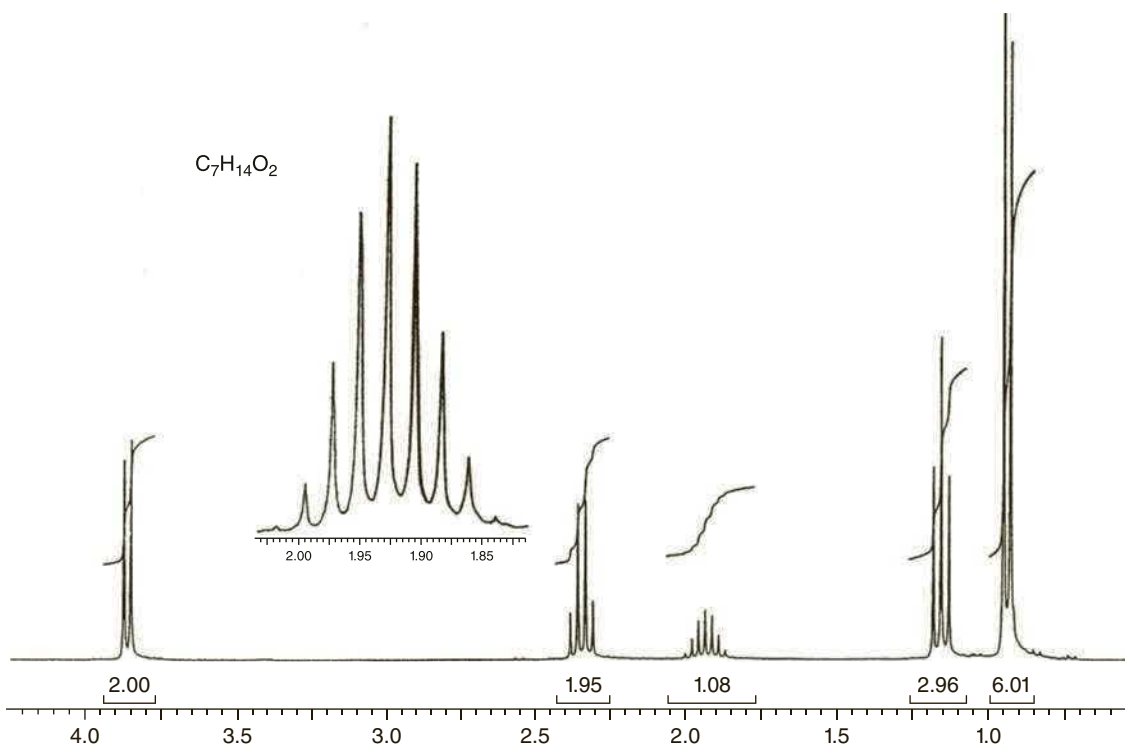

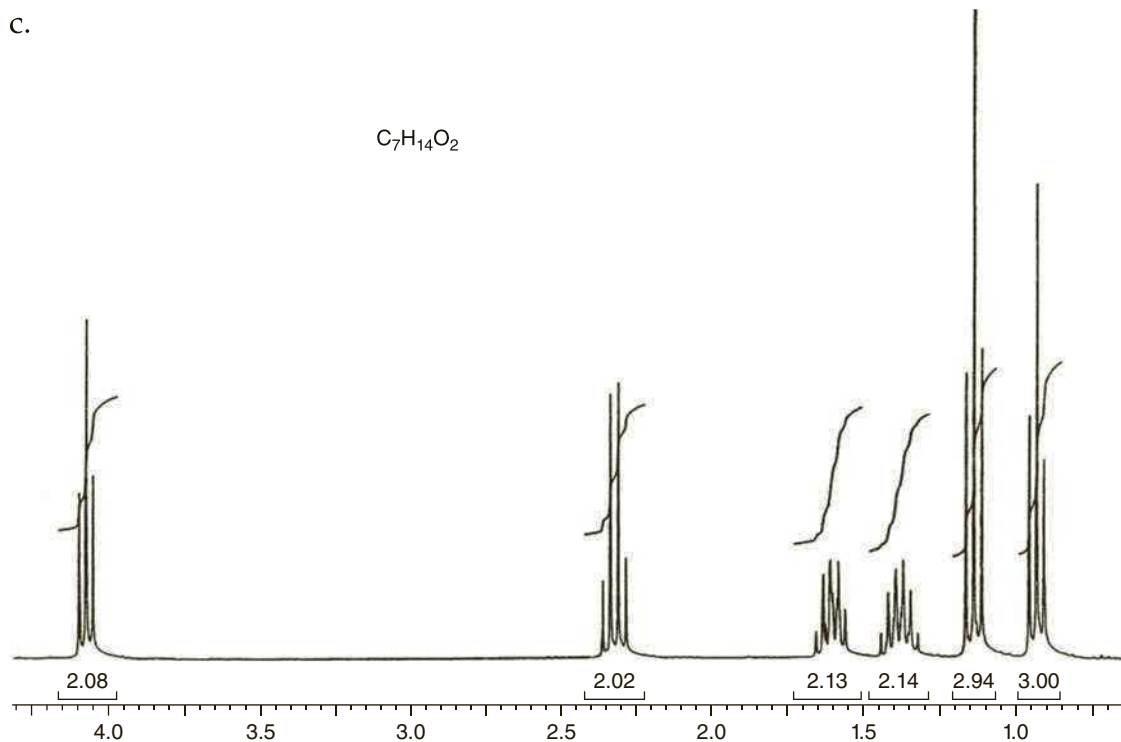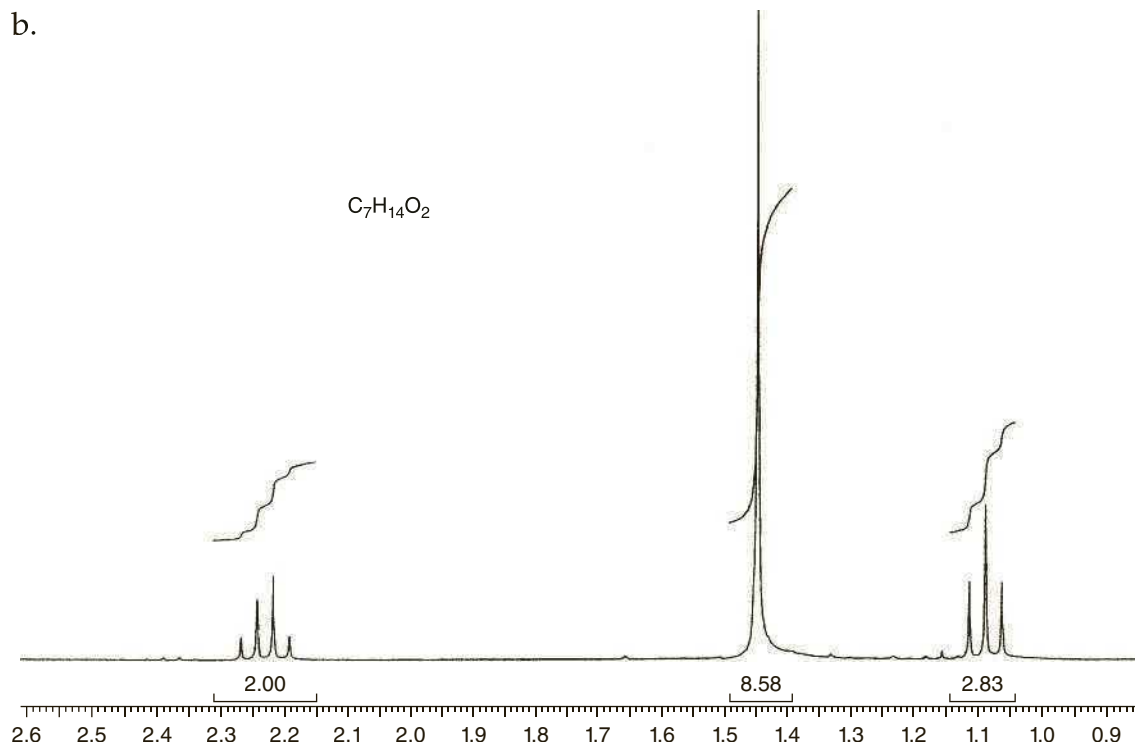


9. The two isomeric carboxylic acids that give the following NMR spectra both have the formula $C_3H_5ClO_2$. Draw their structures.

- a. The broad singlet integrating for one proton that is shown as an inset on the spectrum appears downfield at 11.5 ppm.

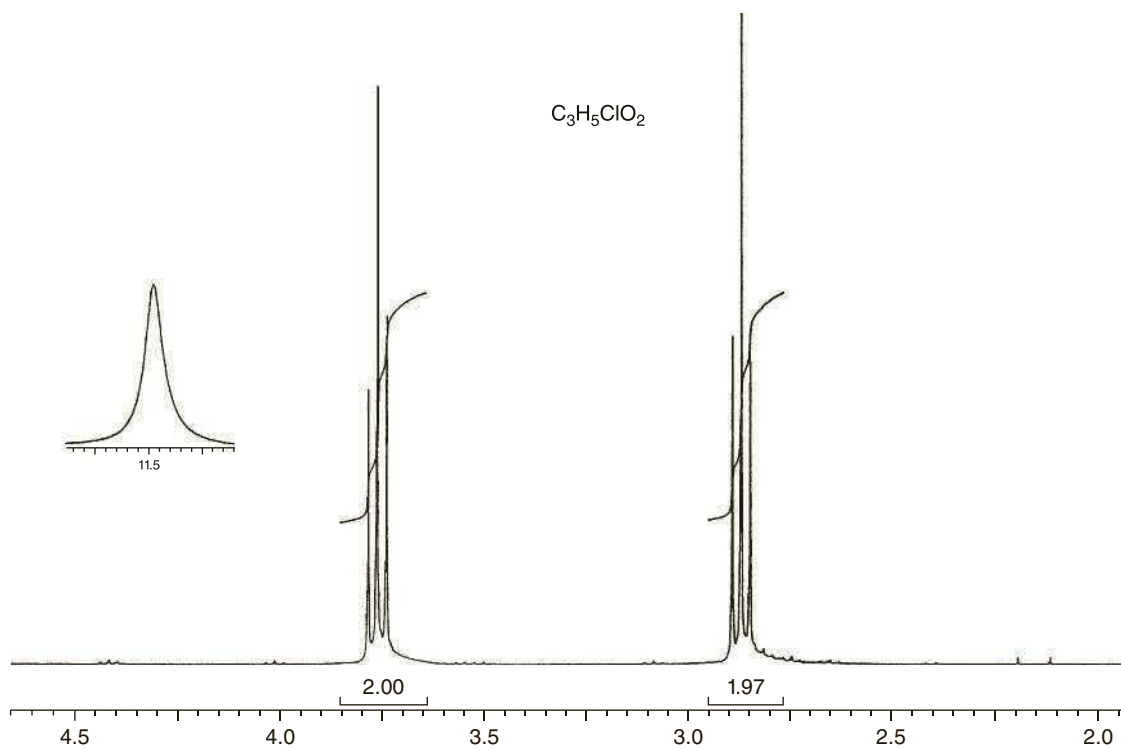


- b. The singlet integrating for one proton that is shown as an inset on the spectrum appears downfield at 12.0 ppm.

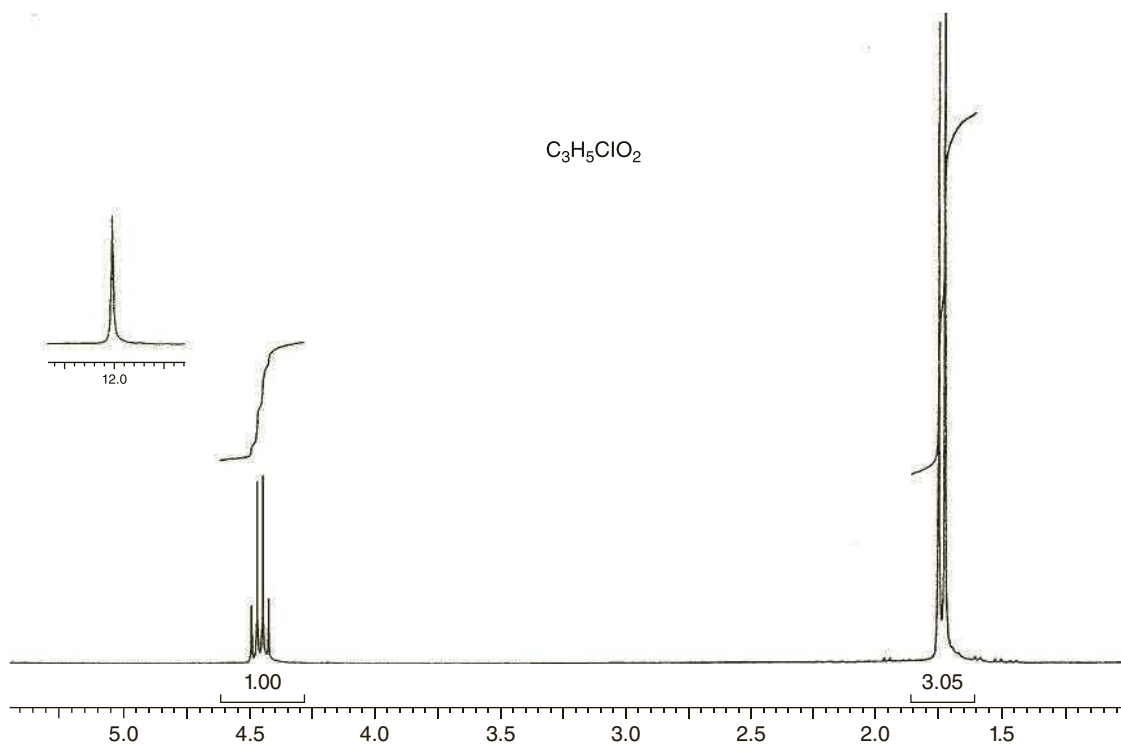


10. The following compounds are isomers with formula $C_{10}H_{12}O$. Their infrared spectra show strong bands near 1715 cm^{-1} and in the range from 1600 cm^{-1} to 1450 cm^{-1} . Draw their structures.

a.

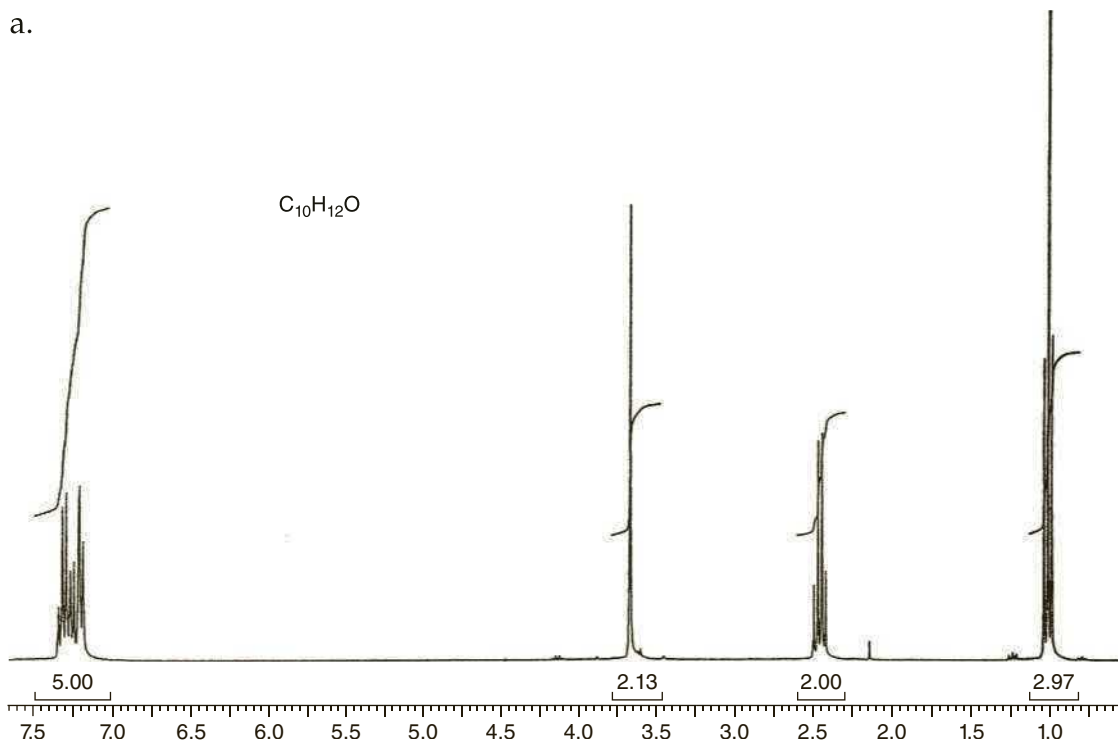

b.

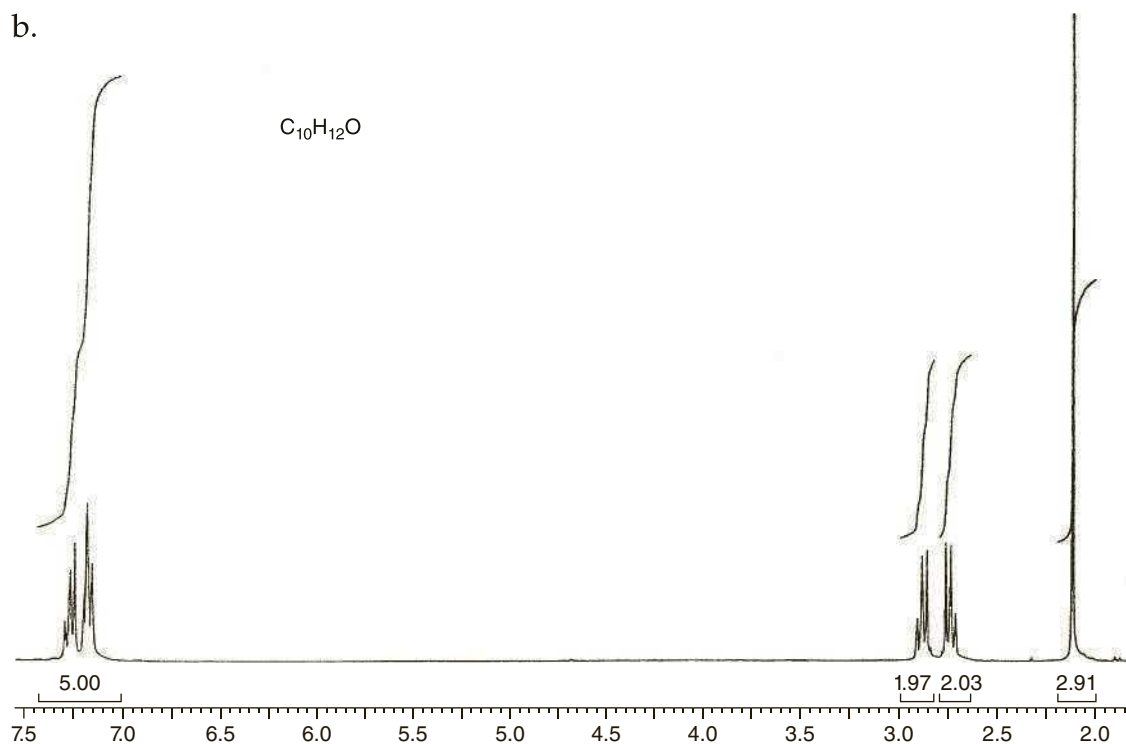

Carbon-13 Nuclear Magnetic Resonance Spectroscopy

Carbon-12, the most abundant isotope of carbon, does not possess spin ($I = 0$); it has both an even atomic number and an even atomic weight. The second principal isotope of carbon, ^{13}C , however, does have the nuclear spin property ($I = \frac{1}{2}$). ^{13}C atom resonances are not easy to observe, due to a combination of two factors. First, the natural abundance of ^{13}C is low; only 1.08% of all carbon atoms are ^{13}C . Second, the magnetic moment μ of ^{13}C is low. For these two reasons, the resonances of ^{13}C are about 6000 times weaker than those of hydrogen. With special Fourier transform (FT) instrumental techniques, which are not discussed here, it is possible to observe ^{13}C nuclear magnetic resonance (carbon-13) spectra on samples that contain only the natural abundance of ^{13}C .

The most useful parameter derived from carbon-13 spectra is the chemical shift. Integrals are unreliable and are not necessarily related to the relative numbers of ^{13}C atoms present in the sample. Hydrogens that are attached to ^{13}C atoms cause spin-spin splitting, but spin-spin interaction between adjacent carbon atoms is rare. With the low natural abundance of carbon-13 (0.0108), the probability of finding two ^{13}C atoms adjacent to one another is extremely low.

Carbon spectra can be used to determine the number of non-equivalent carbons and to identify the types of carbon atoms (methyl, methylene, aromatic, carbonyl, and so on) that may be present in a compound. Thus, carbon NMR provides direct information about the carbon skeleton of a molecule. Because of the low natural abundance of carbon-13 in a sample, it is often necessary to acquire multiple scans over what is needed for proton NMR.

For a given magnetic field strength, the resonance frequency of a ^{13}C nucleus is about one-fourth the frequency required to observe proton resonances. For example, in a 7.05-tesla applied magnetic field, protons are observed at 300 MHz, and ^{13}C nuclei are observed at about 75 MHz.

27.1 Preparing a Sample for Carbon-13 NMR

Technique 26, Section 26.1, describes the technique for preparing samples for proton NMR. Much of what is described there also applies to carbon NMR. There are some differences, however, in determining a carbon spectrum. Fourier transform instruments require a deuterium signal to stabilize (lock) the field. Therefore, the solvents must contain deuterium. Deuterated chloroform, CDCl_3 , is used most commonly for this purpose because of its relatively low cost. Other deuterated solvents may also be used.

Modern FT-NMR spectrometers allow chemists to obtain both the proton and carbon NMR spectra of the same sample in the same NMR tube. After changing several parameters in the program operating the spectrometer, you can obtain both spectra without removing the sample from the probe. The only real difference

is that a proton spectrum may be obtained after a few scans, whereas the carbon spectrum may require 10–100 times more scans.

Tetramethylsilane (TMS) may be added as an internal reference standard, where the chemical shift of the methyl carbon is defined as 0.00 ppm. Alternatively, you may use the center peak of the CDCl_3 pattern, which is found at 77.0 ppm. This pattern can be observed as a small “triplet” near 77.0 ppm in a number of the spectra given in this chapter.

27.2 Carbon-13 Chemical Shifts

An important parameter derived from carbon-13 spectra is the chemical shift. The correlation chart in Figure 27.1 shows typical ^{13}C chemical shifts, listed in parts per million (ppm) from TMS, where the carbons of the methyl groups of TMS (not the hydrogens) are used for reference. Notice that the chemical shifts appear over a range (0–220 ppm) much larger than that observed for protons (0–12 ppm). Because of the very large range of values, nearly every non-equivalent carbon atom in an organic molecule gives rise to a peak with a different chemical shift. Peaks rarely overlap as they often do in proton NMR.

The correlation chart is divided into four sections. Saturated carbon atoms appear at the highest field, nearest to TMS (8–60 ppm). The next section of the chart demonstrates the effect of electronegative atoms (40–80 ppm). The third section includes alkene and aromatic-ring carbon atoms (100–175 ppm). Finally, the fourth section contains carbonyl carbons, which appear at the lowest field values (155–220 ppm).

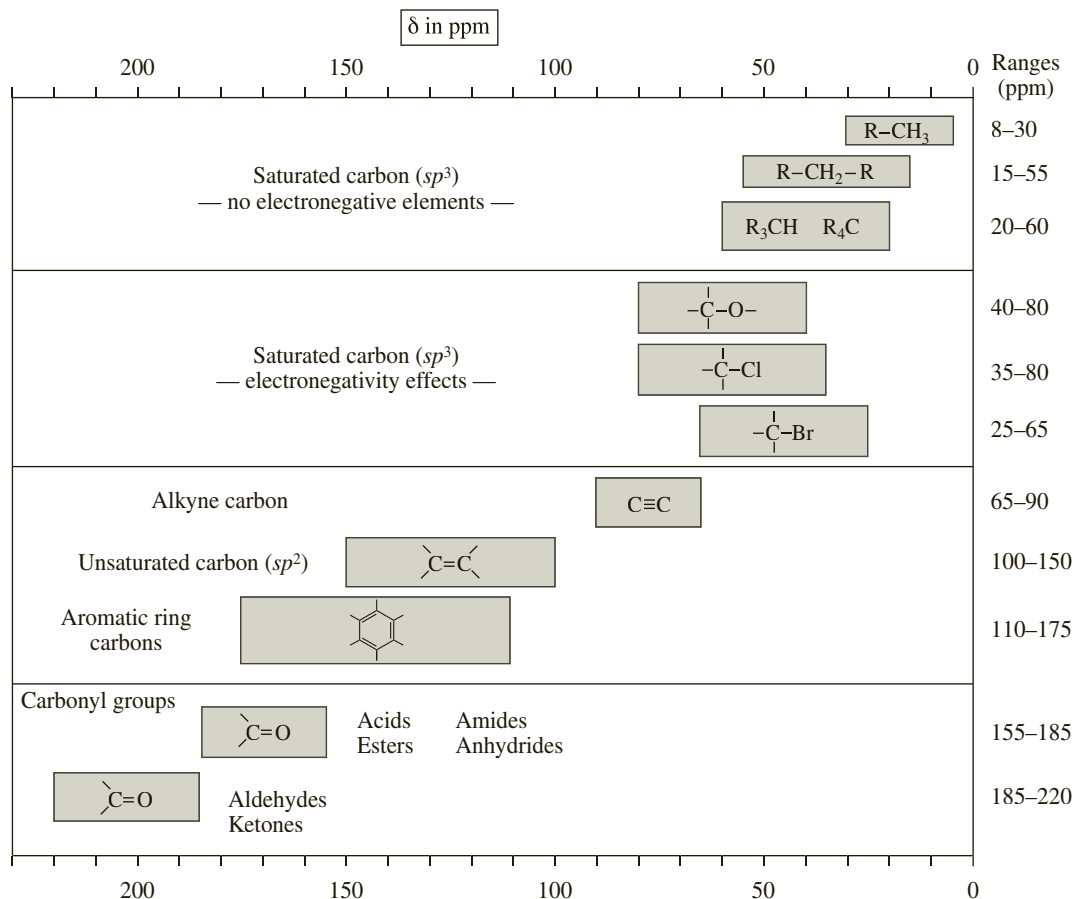
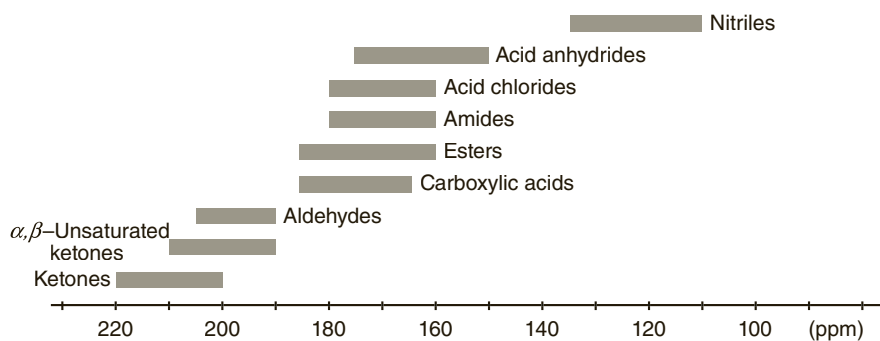


Figure 27.1

A correlation chart for ^{13}C chemical shifts (chemical shifts are listed in parts per million from tetramethylsilane).

**Figure 27.2**

A ^{13}C correlation chart for carbonyl and nitrile functional groups.

Electronegativity, hybridization, and anisotropy all affect ^{13}C chemical shifts in nearly the same fashion as they affect ^1H chemical shifts; however ^{13}C chemical shifts are about 20 times larger. Electronegativity (see Section 26.7) produces the same deshielding effect in carbon NMR as in proton NMR—the electronegative element produces a large downfield shift. The shift is greater for a ^{13}C atom than for a proton because the electronegative atom is directly attached to the ^{13}C atom and the effect occurs through only a single bond, $\text{C}-\text{X}$. With protons, the electronegative atoms are attached to carbon, not hydrogen; the effect occurs through two bonds, $\text{H}-\text{C}-\text{X}$, rather than one.


Analogous with ^1H shifts, changes in hybridization also produce larger shifts for the carbon-13 that is *directly involved* (no bonds) than they do for the hydrogens attached to that carbon (one bond). In ^{13}C NMR, the carbons of carbonyl groups have the largest chemical shifts, due both to sp^2 hybridization and to the fact that an electronegative oxygen is directly attached to the carbonyl carbon, deshielding it even further. Anisotropy (see Section 26.8) is responsible for the large chemical shifts of the carbons in aromatic rings and alkenes.

Notice that the range of chemical shifts is larger for carbon atoms than for hydrogen atoms. Because the factors affecting carbon shifts operate either through one bond or directly on carbon, they are greater than those for hydrogen, which operate through more bonds. As a result, the entire range of chemical shifts becomes larger for ^{13}C (0–220 ppm) than for ^1H (0–12 ppm).

Many of the important functional groups of organic chemistry contain a carbonyl group. In determining the structure of a compound containing a carbonyl group, it is frequently helpful to have some idea of the type of carbonyl group in the unknown. Figure 27.2 illustrates the typical ranges of ^{13}C chemical shifts for some carbonyl-containing functional groups. Although there is some overlap in the ranges, ketones and aldehydes are easy to distinguish from the other types. Chemical shift data for carbonyl carbons are particularly powerful when combined with data from an infrared spectrum.

27.3 Proton-Coupled ^{13}C Spectra—Spin-Spin Splitting of Carbon-13 Signals

Unless a molecule is artificially enriched by synthesis, the probability of finding two ^{13}C atoms in the same molecule is low. The probability of finding two ^{13}C atoms adjacent to each other in the same molecule is even lower. Therefore, we rarely observe homonuclear (carbon–carbon) spin–spin splitting patterns where the interaction occurs between two ^{13}C atoms. However, the spins of protons attached directly to ^{13}C atoms do interact with the spin of carbon and cause the carbon signal to be split according to the $n + 1$ rule. This is **heteronuclear** (carbon–hydrogen) coupling involving two different types of atoms. With ^{13}C NMR, we generally examine splitting

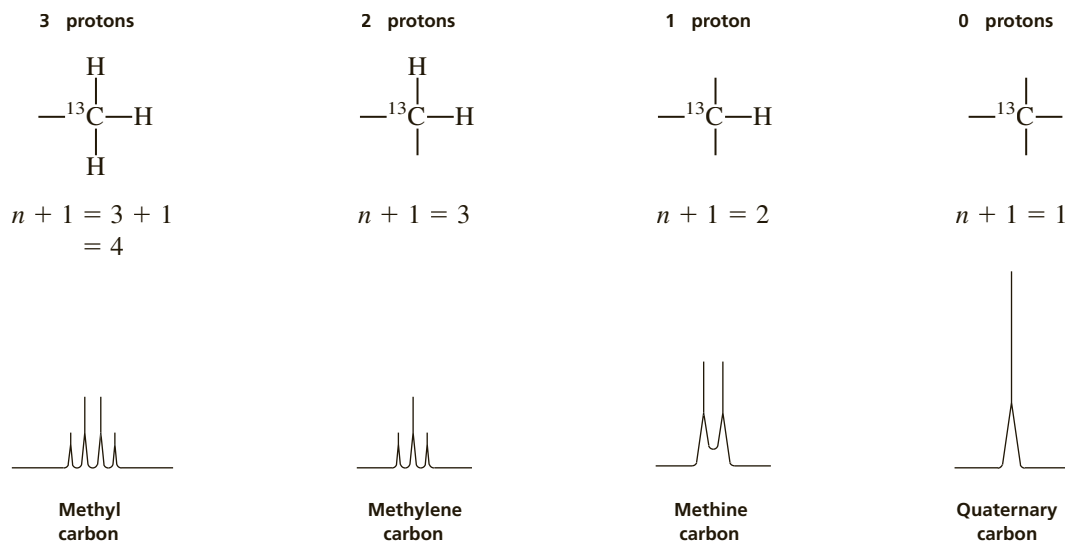**Figure 27.3**

The effect of attached protons on ^{13}C resonances.

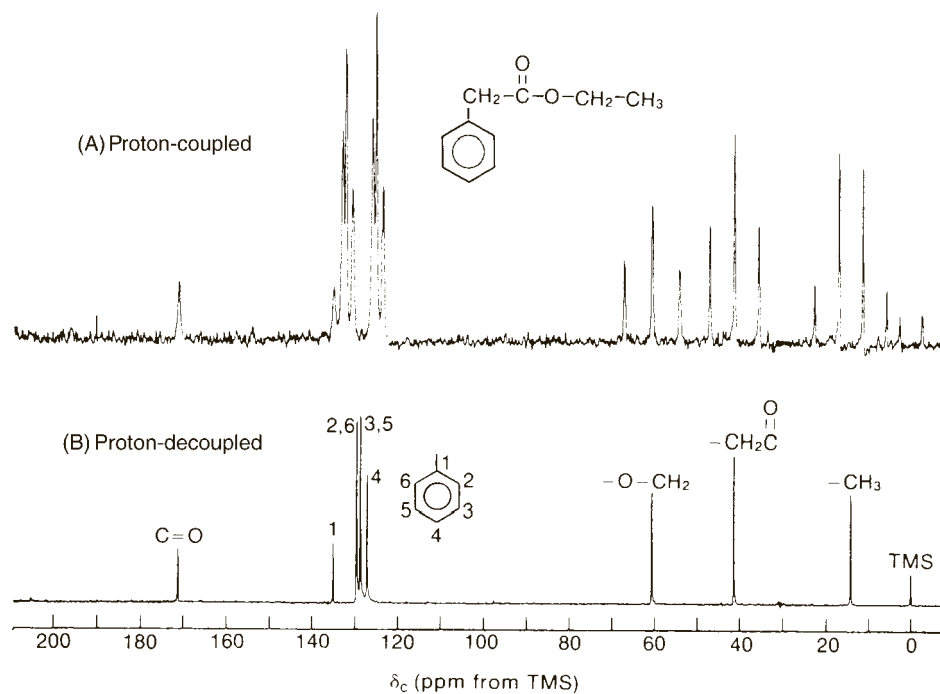
that arises from the protons *directly attached* to the carbon atom being studied. This is a one-bond coupling. In proton NMR, the most common splittings are *homonuclear* (hydrogen–hydrogen), which occur between protons attached to *adjacent* carbon atoms. In these cases, the interaction is a three-bond coupling, $\text{H}-\text{C}-\text{C}-\text{H}$.

Figure 27.3 illustrates the effect of protons directly attached to a ^{13}C atom. The $n + 1$ rule predicts the degree of splitting in each case. The resonance of a ^{13}C atom with three attached protons, for instance, is split into a quartet ($n + 1 = 3 + 1 = 4$). Because the hydrogens are directly attached to the carbon-13 (one-bond couplings), the coupling constants for this interaction are quite large, with J values of about 100 Hz to 250 Hz. Compare the typical three-bond $\text{H}-\text{C}-\text{C}-\text{H}$ couplings that are common in NMR spectra, which have J values of about 4 Hz to 18 Hz.

It is important to note while examining Figure 27.3 that you are not “seeing” protons directly when looking at a ^{13}C spectrum (proton resonances occur at frequencies outside the range used to obtain ^{13}C spectra); you are observing only the effect of the protons on ^{13}C atoms. Also remember that we cannot observe ^{12}C , because it is NMR inactive.

Spectra that show the spin–spin splitting, or coupling, between carbon-13 and the protons directly attached to it are called **proton-coupled spectra**. Figure 27.4 A is the proton-coupled ^{13}C NMR spectrum of ethyl phenylacetate. In this spectrum, the first quartet downfield from TMS (14.2 ppm) corresponds to the carbon of the methyl group. It is split into a quartet ($J = 127$ Hz) by the three attached hydrogen atoms ($^{13}\text{C}-\text{H}$, one-bond couplings). In addition, although it cannot be seen on the scale of this spectrum (an expansion must be used), each of the quartet lines is split into a closely spaced triplet ($J = \text{ca. } 1$ Hz). This additional fine splitting is caused by the two protons on the adjacent $-\text{CH}_2-$ group. These are two-bond couplings ($\text{H}-\text{C}-^{13}\text{C}$) of a type that occurs commonly in ^{13}C spectra, with coupling constants that are generally quite small ($J = 0\text{--}2$ Hz) for systems with carbon atoms in an aliphatic chain. Because of their small size, these couplings are frequently ignored in the routine analysis of spectra, with greater attention being given to the larger one-bond splittings seen in the quartet itself.

There are two $-\text{CH}_2-$ groups in ethyl phenylacetate. The one corresponding to the ethyl $-\text{CH}_2-$ group is found farther downfield (60.6 ppm), as this carbon is deshielded by the attached oxygen. It is a triplet because of the two attached

**Figure 27.4**

Ethyl phenylacetate. (A) The proton-coupled ^{13}C NMR spectrum (20 MHz). (B) The proton-decoupled ^{13}C spectrum (20 MHz). (From Moore, J. A., Dalrymple, D. L., and Rodig, O. R. *Experimental Methods in Organic Chemistry*, 3rd ed. [W. B. Saunders: Philadelphia, 1982])

hydrogens (one-bond couplings). Again, although it is not seen in this unexpanded spectrum, the three hydrogens on the adjacent methyl group finely split each of the triplet peaks into a quartet. The benzyl $-\text{CH}_2-$ carbon is the intermediate triplet (41.4 ppm). Farthest downfield is the carbonyl-group carbon (171.1 ppm). On the scale of this presentation, it is a singlet (no directly attached hydrogens), but because of the adjacent benzyl $-\text{CH}_2-$ group, it is actually split finely into a triplet. The aromatic ring carbons also appear in the spectrum, and they have resonances in the range from 127 ppm to 136 ppm. Section 27.7 will discuss aromatic ring ^{13}C resonances.

Proton-coupled spectra for large molecules are often difficult to interpret. The multiplets from different carbons commonly overlap because the $^{13}\text{C}-\text{H}$ coupling constants are frequently larger than the chemical shift differences of the carbons in the spectrum. Sometimes, even simple molecules such as ethyl phenylacetate (Figure 27.4A) are difficult to interpret. Proton decoupling, which is discussed in the next section, avoids this problem.

27.4 Proton-Decoupled ^{13}C Spectra

By far, the great majority of ^{13}C NMR spectra are obtained as proton-decoupled spectra. The decoupling technique obliterates all interactions between protons and ^{13}C nuclei; therefore, only singlets are observed in a decoupled ^{13}C NMR spectrum. Although this technique simplifies the spectrum and avoids overlapping multiplets, it has the disadvantage that the information on attached hydrogens is lost.

Proton **decoupling** is accomplished in the process of determining a ^{13}C NMR spectrum by simultaneously irradiating all of the protons in the molecule with a broad spectrum of frequencies in the proper range for protons. Modern NMR

spectrometers provide a second, tunable radio-frequency generator, the decoupler, for this purpose. Irradiation causes the protons to become saturated, and they undergo rapid upward and downward transitions, among all their possible spin states. These rapid transitions decouple any spin–spin interactions between the hydrogens and the ^{13}C nuclei being observed. In effect, all spin interactions are averaged to zero by the rapid changes. The carbon nucleus “senses” only one average spin state for the attached hydrogens rather than two or more distinct spin states.

Figure 27.4B is a proton-decoupled spectrum of ethyl phenylacetate. The proton-coupled spectrum (see Figure 27.4A) was discussed in Section 27.3. It is interesting to compare the two spectra to see how the proton-decoupling technique simplifies the spectrum. Every chemically and magnetically distinct carbon gives only a single peak. Notice, however, that the two *ortho* ring carbons (carbons 2 and 6) and the two *meta* ring carbons (carbons 3 and 5) are equivalent by symmetry and that each pair gives only a single peak.

Figure 27.5 is a second example of a proton-decoupled spectrum. Notice that the spectrum shows three peaks corresponding to the exact number of carbon atoms in 1-propanol. If there are no equivalent carbon atoms in a molecule, a ^{13}C peak will be observed for *each* carbon. Notice also that the assignments given in Figure 27.5 are consistent with the values in the chemical shift chart (see Figure 27.1). The carbon atom closest to the electronegative oxygen is farthest downfield, and the methyl carbon is at highest field.

The three-peak pattern centered at $\delta = 77$ ppm is due to the solvent CDCl_3 . This pattern results from the coupling of a deuterium (^2H) nucleus to the ^{13}C nucleus. Often, the CDCl_3 pattern is used as an internal reference in place of TMS.


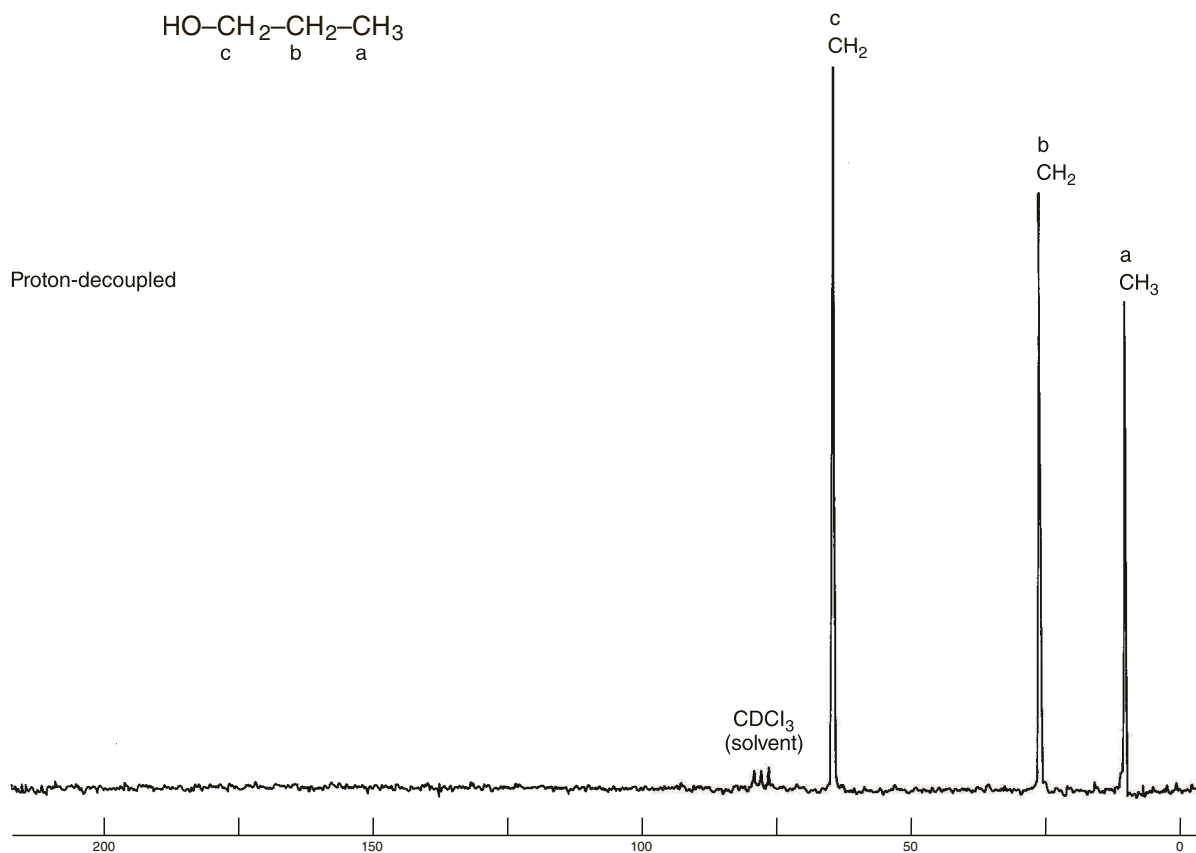


Figure 27.5

The proton-decoupled ^{13}C NMR spectrum of 1-propanol (22.5 MHz).

27.5 Some Sample Spectra—Equivalent Carbons

Equivalent ^{13}C atoms appear at the same chemical shift value. Figure 27.6 shows the proton-decoupled carbon spectrum for 2,2-dimethylbutane. The three methyl groups at the left side of the molecule are equivalent by symmetry.

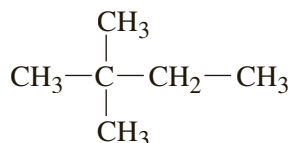

Although this compound has a total of six carbons, there are only four peaks in the ^{13}C NMR spectrum. The ^{13}C atoms that are equivalent appear at the same chemical shift. The single methyl carbon, a, appears at highest field (9 ppm), and the three equivalent methyl carbons, b, appear at 29 ppm. The quaternary carbon, c, gives rise to the small peak at 30 ppm, and the methylene carbon, d, appears at 37 ppm. The relative sizes of the peaks are related, in part, to the number of each type of carbon atom present in the molecule. For example, notice in Figure 27.6 that the peak at 29 ppm (b) is much larger than the others. This peak is generated by three carbons. The quaternary carbon at 30 ppm (c) is very weak. Because no hydrogens are attached to this carbon, there is very little nuclear Overhauser enhancement (NOE) (see Section 27.6). Without attached hydrogen atoms, relaxation times are also longer than for other carbon atoms. Quaternary carbons, those with no hydrogens attached, frequently appear as weak peaks in proton-decoupled ^{13}C NMR spectra (see Section 27.6).

Figure 27.7 is a proton-decoupled ^{13}C spectrum of cyclohexanol. This compound has a plane of symmetry passing through its hydroxyl group, and it shows only four carbon resonances. Carbons a and c are doubled due to symmetry and give rise to larger peaks than carbons b and d. Carbon d, bearing the hydroxyl group, is deshielded by oxygen and has its peak at 70.0 ppm. Notice that this peak has the lowest intensity of all of the peaks. Its intensity is lower than that of carbon b in part because the carbon d peak receives the least amount of NOE; there is only one hydrogen attached to the hydroxyl carbon, whereas each of the other carbons has two hydrogens.

A carbon attached to a double bond is deshielded due to its sp^2 hybridization and some diamagnetic anisotropy. This effect can be seen in the ^{13}C NMR spectrum of cyclohexene (Figure 27.8). Cyclohexene has a plane of symmetry that runs

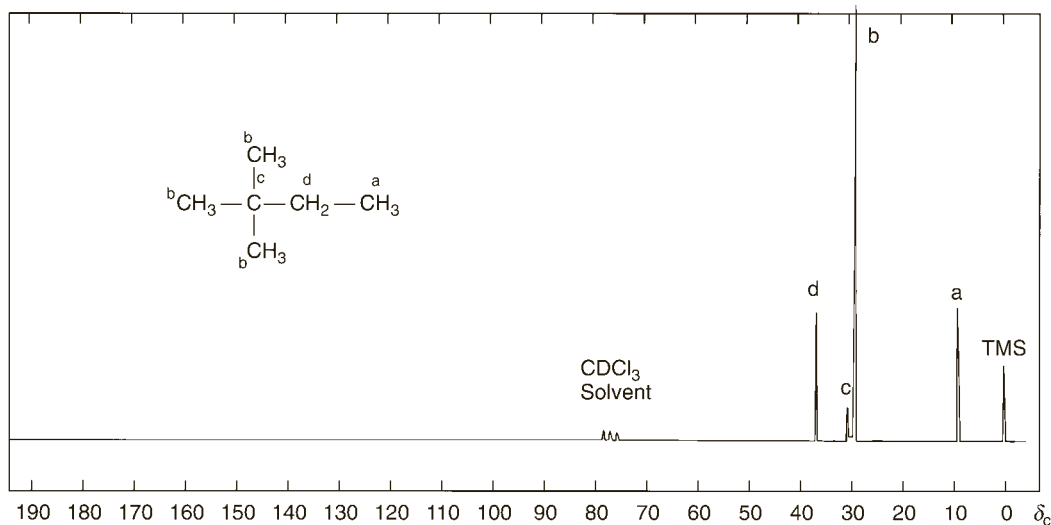


Figure 27.6
The proton-decoupled ^{13}C NMR spectrum of 2,2-dimethylbutane.

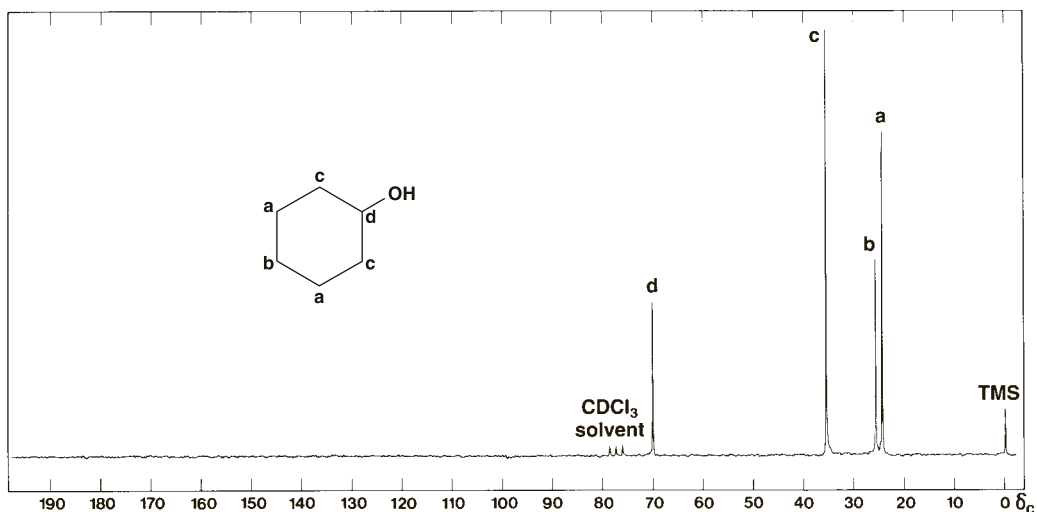


Figure 27.7
The proton-decoupled ^{13}C NMR spectrum of cyclohexanol.

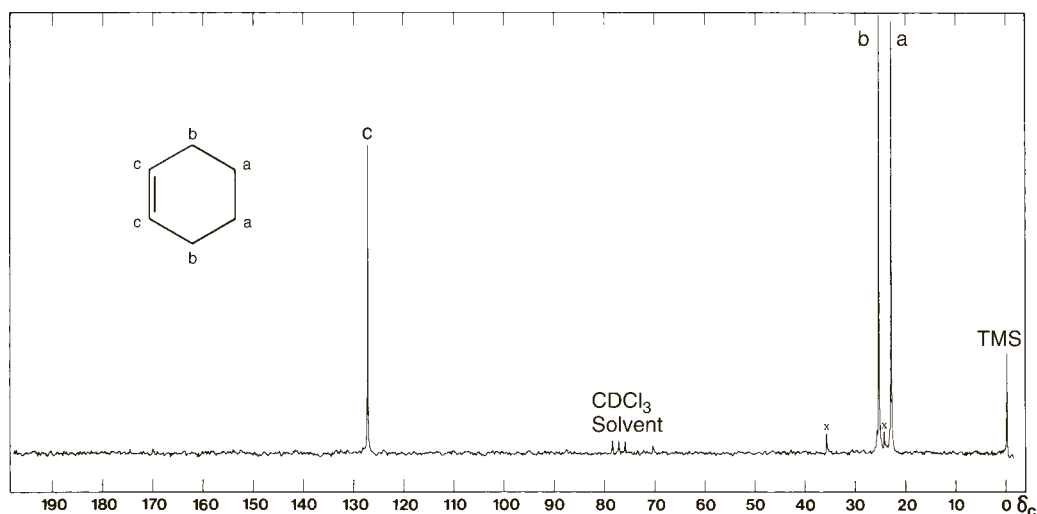


Figure 27.8
The proton-decoupled ^{13}C NMR spectrum of cyclohexene. (The peaks marked with an x results from impurities.)

perpendicular to the double bond. As a result, we observe only three absorption peaks. There are two of each type of sp^3 carbon. Each of the double-bond carbons *c* has only one hydrogen, whereas each of the remaining carbons has two. As a result of a reduced NOE, the double-bond carbons (127 ppm) have a lower-intensity peak in the spectrum.

In Figure 27.9, the spectrum of cyclohexanone, the carbonyl carbon has the lowest intensity. This is due not only to reduced NOE (no hydrogen attached) but also to the long relaxation time of the carbonyl carbon (see Section 27.6). Notice also that Figure 27.2 predicts the large chemical shift for this carbonyl carbon (211 ppm).

27.6 Nuclear Overhauser Enhancement (NOE)

When we obtain a proton-decoupled ^{13}C spectrum, the intensities of many of the carbon resonances increase significantly above those observed in a proton-coupled experiment. Carbon atoms with hydrogen atoms directly attached are enhanced the most, and the enhancement increases (but not always linearly) as more hydrogens


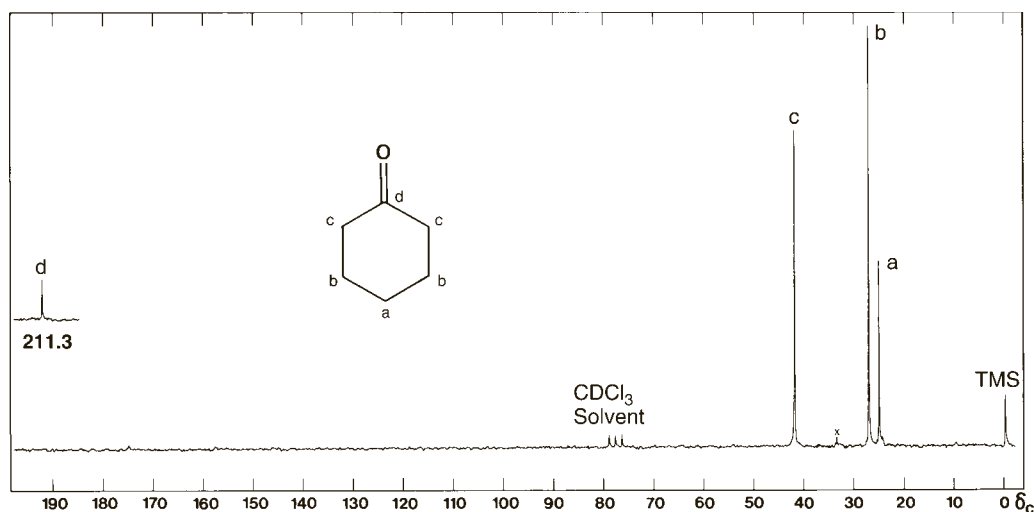


Figure 27.9

The proton-decoupled ^{13}C NMR spectrum of cyclohexanone. (The peak marked with an x is an impurity.)

are attached. This effect is known as the nuclear Overhauser effect, and the degree of increase in the signal is called the **nuclear Overhauser enhancement (NOE)**. Thus, we expect that the intensity of the carbon peaks should increase in the following order in a typical carbon-13 NMR spectrum:

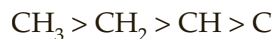


Carbon atom relaxation times influence the intensity of peaks in a spectrum. When more protons are attached to a carbon atom, relaxation times become shorter, resulting in more intense peaks. We expect methyl and methylene groups to be relatively more intense than the intensity observed for quaternary carbon atoms where there are no attached protons. Thus, a weak-intensity peak is observed for the quaternary carbon atom at 30 ppm in 2,2-dimethylbutane (see Figure 27.6).

In addition, weak carbonyl carbon peaks are observed at 171 ppm in ethyl phenylacetate (see Figure 27.4) and at 211 ppm in cyclohexanone (see Figure 27.9).

27.7 Compounds with Aromatic Rings

Compounds with carbon-carbon double bonds or aromatic rings give rise to chemical shifts from 100 ppm to 175 ppm. Because relatively few other peaks appear in this range, a great deal of useful information is available when peaks appear here.

A **monosubstituted** benzene ring shows *four* peaks in the aromatic carbon area of a proton-decoupled ^{13}C spectrum, because the *ortho* and *meta* carbons are doubled by symmetry. Often the carbon with no protons attached, the *ipso* carbon, has a very weak peak due to a long relaxation time and a weak NOE. In addition, there are two larger peaks for the doubled *ortho* and *meta* carbons and a medium-sized peak for the *para* carbon. In many cases, it is not important to be able to assign all of the peaks precisely. In the example of toluene, shown in Figure 27.10, notice that carbons c and d are not easy to assign by inspection of the spectrum.

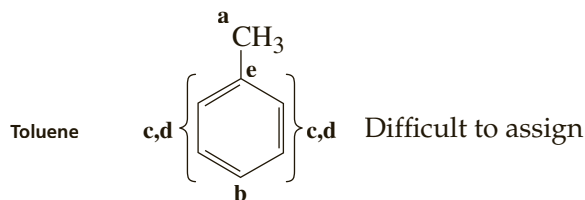

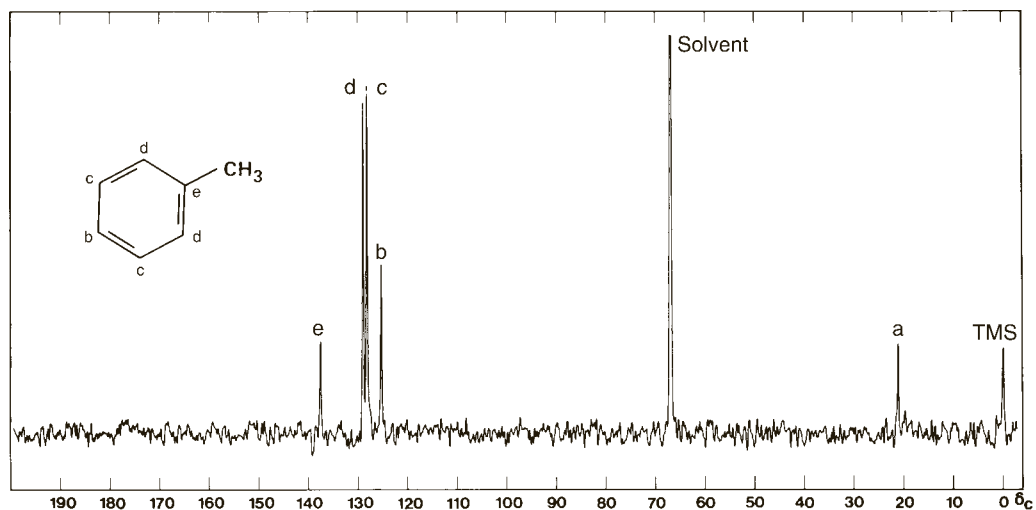


Figure 27.10

The proton-decoupled ^{13}C NMR spectrum of toluene.

In a proton-coupled ^{13}C spectrum, a monosubstituted benzene ring shows three doublets and one singlet. The singlet arises from the *ipso* carbon, which has no attached hydrogen. Each of the other carbons in the ring (*ortho*, *meta*, and *para*) has one attached hydrogen and yields a doublet.

Figure 27.4B is the proton-decoupled spectrum of ethyl phenylacetate, with the assignments noted next to the peaks. Notice that the aromatic ring region shows four peaks between 125 ppm and 135 ppm, consistent with a monosubstituted ring. There is one peak for the methyl carbon (13 ppm) and there are two peaks for the methylene carbons. One of the methylene carbons is directly attached to an electronegative oxygen atom and appears at 61 ppm, and the other is more shielded (41 ppm). The carbonyl carbon (an ester) has resonance at 171 ppm. All of the carbon chemical shifts agree with the values in the correlation chart (Figure 27.1).

Depending on the mode of substitution, a symmetrically **disubstituted** benzene ring can show two, three, or four peaks in the proton-decoupled ^{13}C spectrum. The following drawings illustrate this for the isomers of dichlorobenzene.

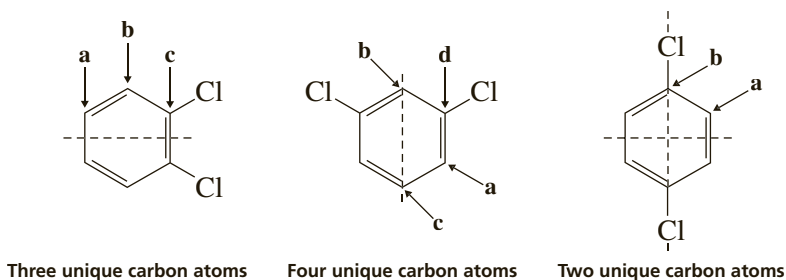


Figure 27.11 shows the spectra of all three dichlorobenzenes, each of which has the number of peaks consistent with the analysis just given. You can see that ^{13}C NMR spectroscopy is very useful in the identification of isomers.

Most other polysubstitution patterns on a benzene ring yield six peaks in the proton-decoupled ^{13}C NMR spectrum, one for each carbon. However, when identical substituents are present, watch carefully for planes of symmetry that may reduce the number of peaks.

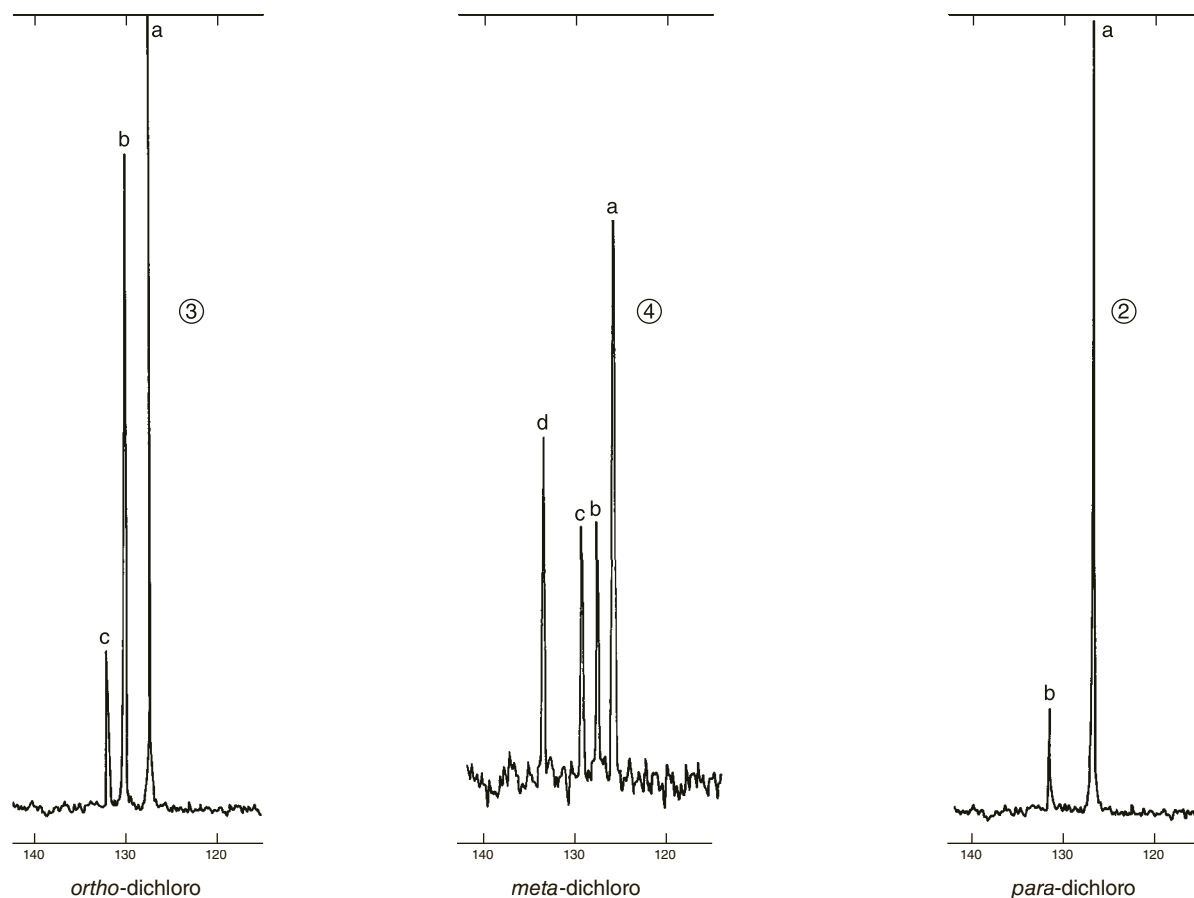


Figure 27.11
The proton-decoupled ^{13}C NMR spectra of the three isomers of dichlorobenzene (25 MHz).

REFERENCES

Textbooks

- Friebolin, H. *Basic One- and Two-Dimensional NMR Spectroscopy*, 5th ed. VCH Publishers: New York, 2010.
- Gunther, H. *NMR Spectroscopy*, 2nd ed. John Wiley & Sons: New York, 1995.
- Levy, G. C. *Topics in Carbon-13 Spectroscopy*. John Wiley & Sons: New York, 1984.
- Levy, G. C.; Lichter, R. L.; Nelson, G. L. *Carbon-13 Nuclear Magnetic Resonance Spectroscopy*, 2nd ed. John Wiley & Sons: New York, 1980.
- Macomber, R. S. *A Complete Introduction to Modern NMR Spectroscopy*. John Wiley & Sons: New York, 2008.
- Macomber, R. S. *NMR Spectroscopy—Essential Theory and Practice*. College Outline Series, Harcourt Brace Jovanovich: New York, 1988.
- Pavia, D. L.; Lampman, G. M.; Kriz, G. S.; and Vyvyan, J. R. *Introduction to Spectroscopy*, 5th ed. Cengage Learning: Stamford, CT, 2015.
- Sanders, J. K. M.; Hunter, B. K. *Modern NMR Spectroscopy—A Guide for Chemists*, 2d ed. Oxford University Press: Oxford, England, 1993.
- Silverstein, R. M.; Webster, F. X.; Kiemle, D.; Bryce, D. L. *Spectrometric Identification of Organic Compounds*, 8th ed. John Wiley & Sons: New York, 2015.

Compilations of Spectra

- Johnson, L. F.; Jankowski, W. C. *Carbon-13 NMR Spectra: A Collection of Assigned, Coded, and Indexed Spectra, 25 MHz*. Wiley-Interscience: New York, 1972.

Pouchert, C. J.; Behnke, J. *The Aldrich Library of ^{13}C and ^1H FT-NMR Spectra, 75 and 300 MHz*. Aldrich Chemical Company: Milwaukee, WI, 1993.

Pretsch, E.; Clerc, T.; Seibl, J.; Simon, W. *Tables of Spectral Data for Structure Determination of Organic Compounds*, 2nd ed. Springer-Verlag: Berlin and New York, 1989. Translated from the German by K. Biemann.

Web Sites

<http://www.aist.go.jp/RIODB/SDBS/menu-e.html>












Integrated Spectral DataBase System for Organic Compounds, National Institute of Materials and Chemical Research, Tsukuba, Ibaraki 305-8565, Japan. This database includes infrared, mass spectra, and NMR data (proton and carbon-13) for a number of compounds.

<http://www.chem.ucla.edu/~webspectra>



UCLA Department of Chemistry and Biochemistry in connection with Cambridge University Isotope Laboratories maintains a Web site, WebSpectra, that provides NMR and IR spectroscopy problems for students to interpret. They provide links to other sites with problems for students to solve.



PROBLEMS

1. Predict the number of peaks that you would expect in the proton-decoupled ^{13}C spectrum of each of the following compounds. Problems 1a and 1b are provided as examples. Dots are used to show the non-equivalent carbon atoms in these two examples.

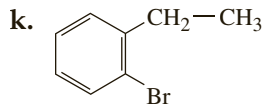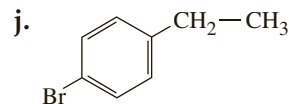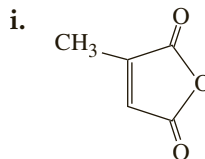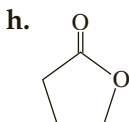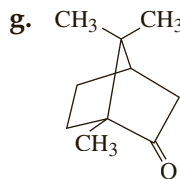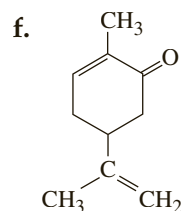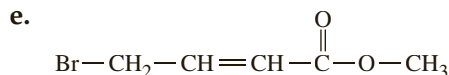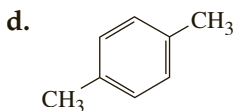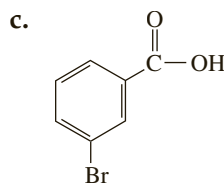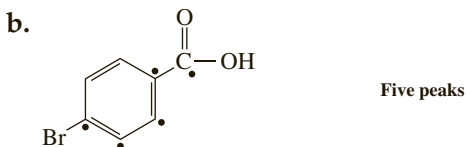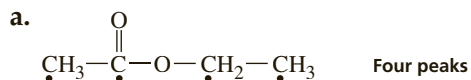

2. Following are the ^1H and ^{13}C spectra for two isomeric bromoalkanes (A and B) with formula $\text{C}_4\text{H}_9\text{Br}$. Integral curves are drawn on the spectra, along with relative integral values provided just above the scale and under each set of peaks. These numbers indicate the relative number of protons assigned to each pattern. Remember that these integral values are approximate. You will need to round the values off to the nearest whole number. Also, in some cases, the lowest whole-number ratios are given. In that case, the values provided may need to be multiplied by two or three in order to obtain the actual number of protons in each pattern.

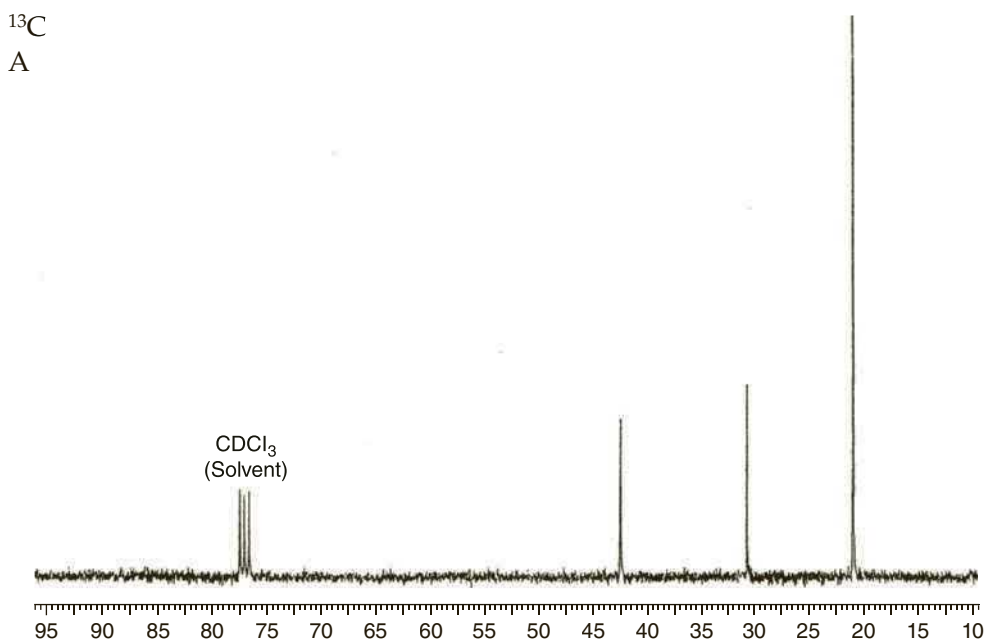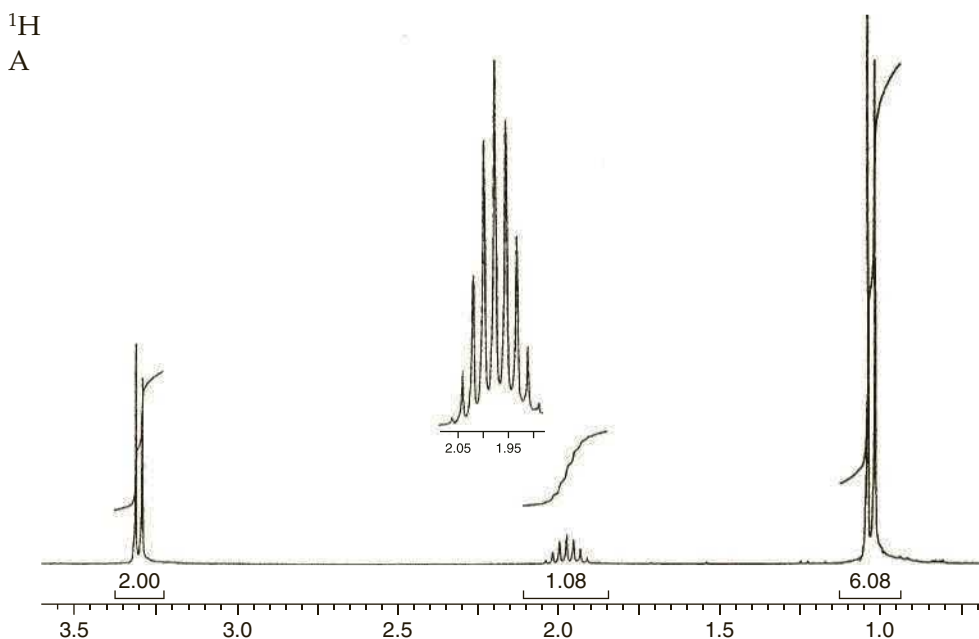


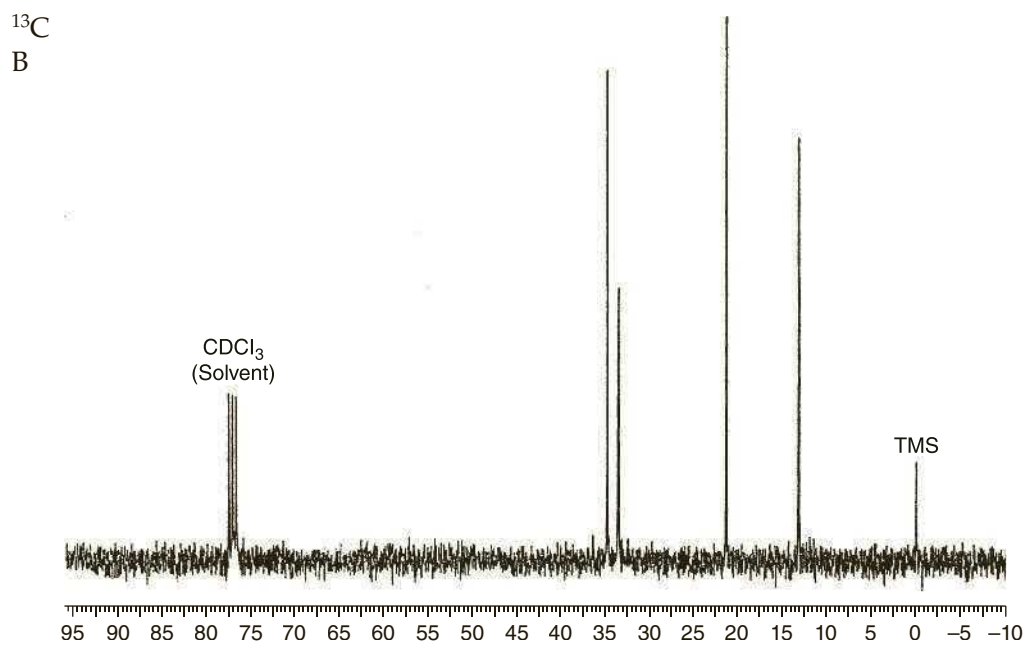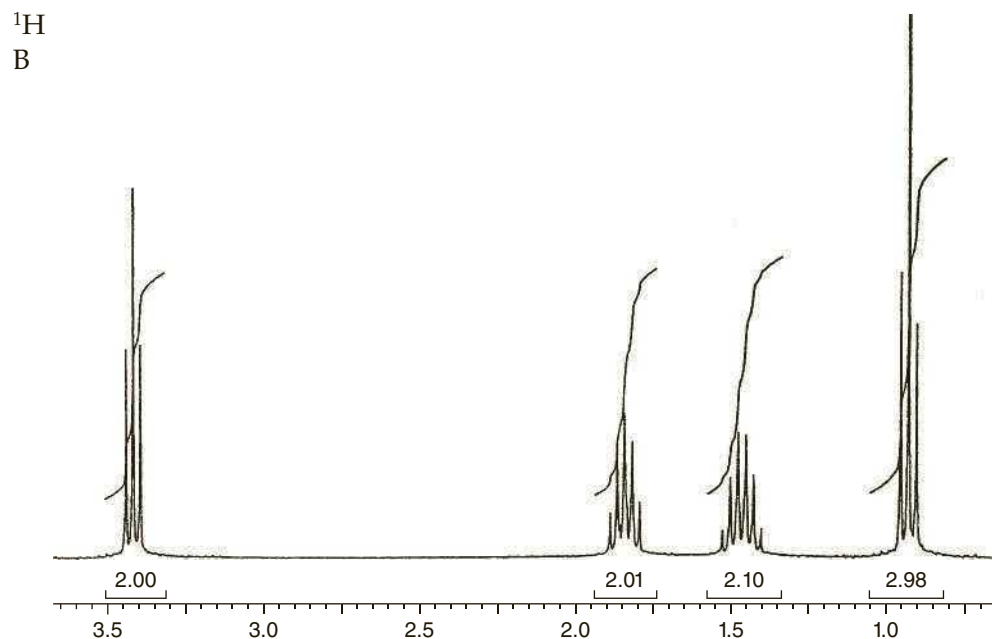


3. Following are the ^1H and ^{13}C spectra for each of three isomeric ketones (A, B, and C) with formula $\text{C}_7\text{H}_{14}\text{O}$. Integral curves are drawn on the spectra, along with relative integral values provided just above the scale and under each set of peaks. These numbers indicate the relative number of protons assigned to each pattern. Remember that these integral values are approximate. You will need to round the values off to the nearest whole number. Also, in some cases, the



lowest whole-number ratios are given. In that case, the values provided may need to be multiplied by two or three in order to obtain the actual number of protons in each pattern.



^1H
A

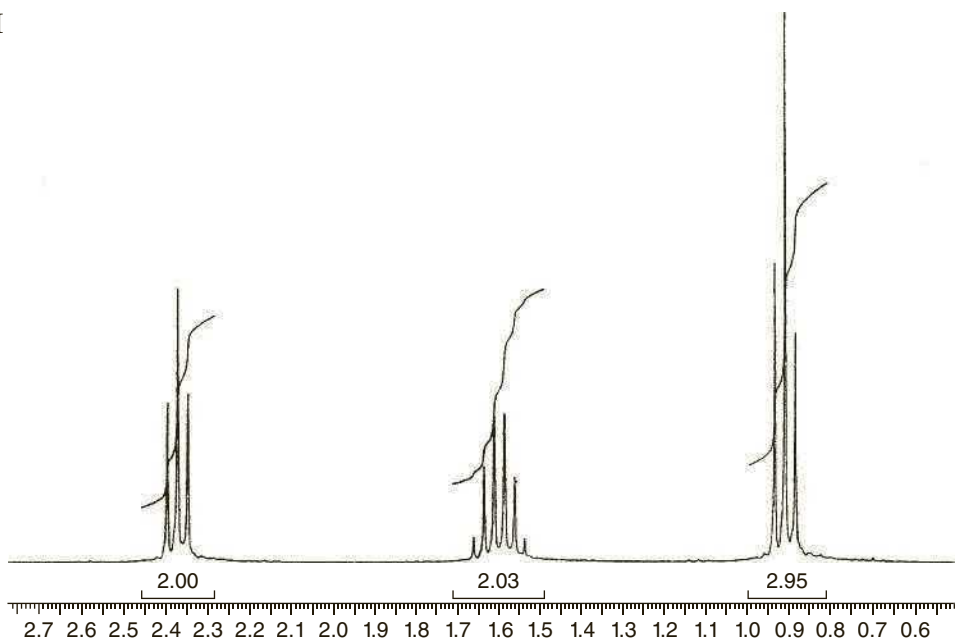

^{13}C
A

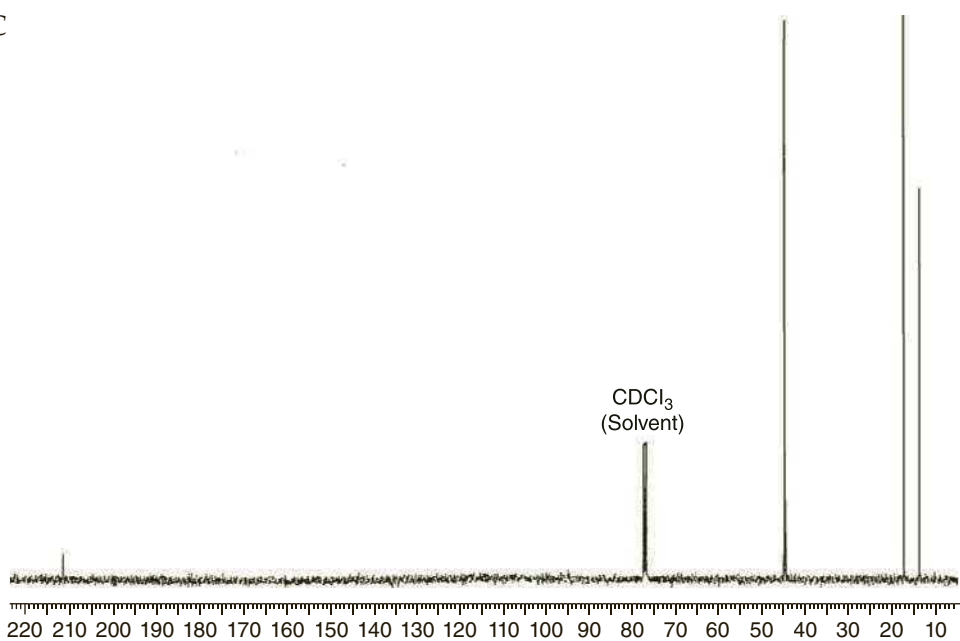

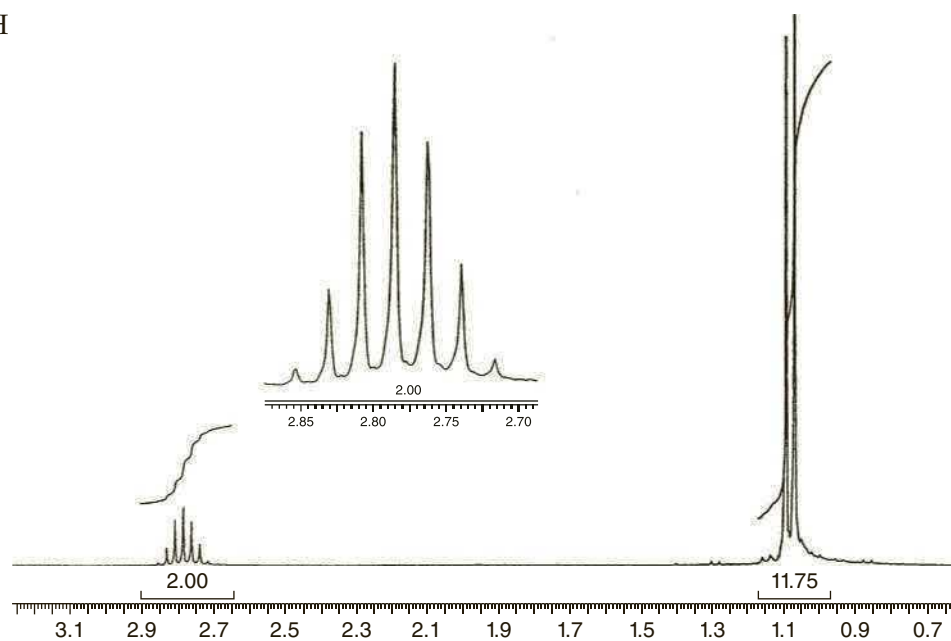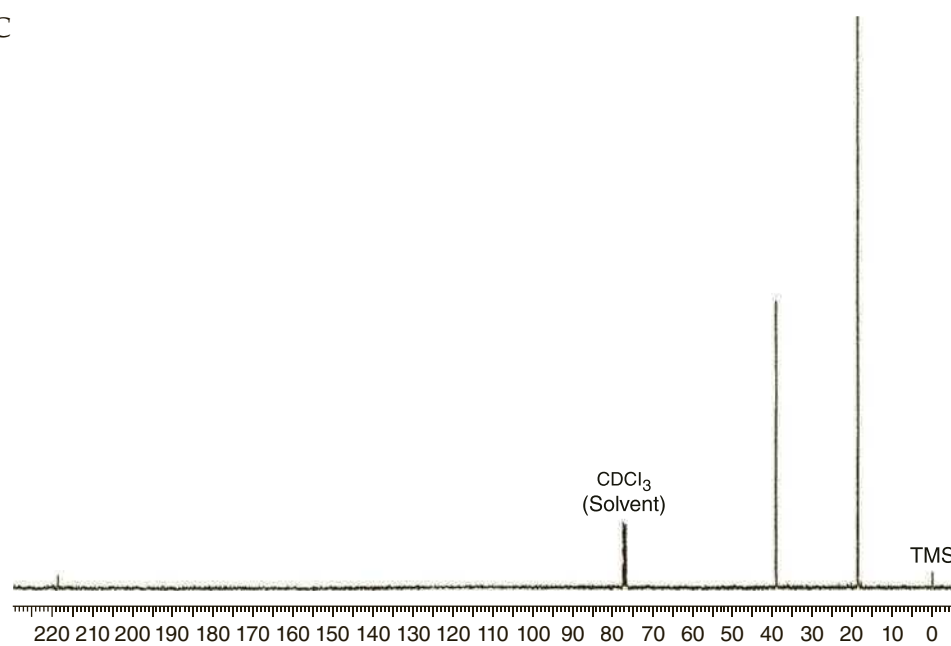^1H
B ^{13}C
B

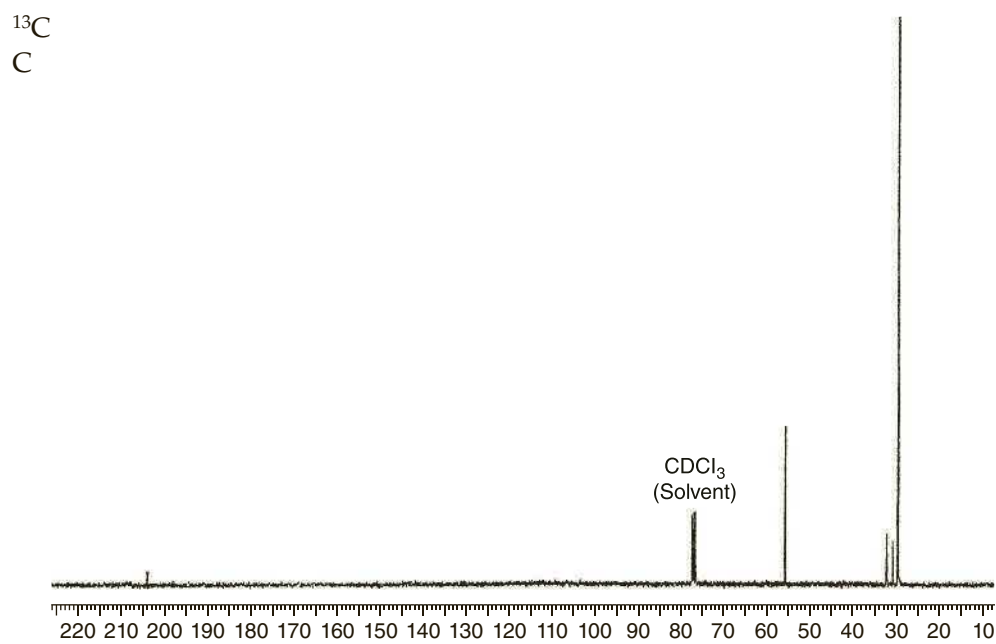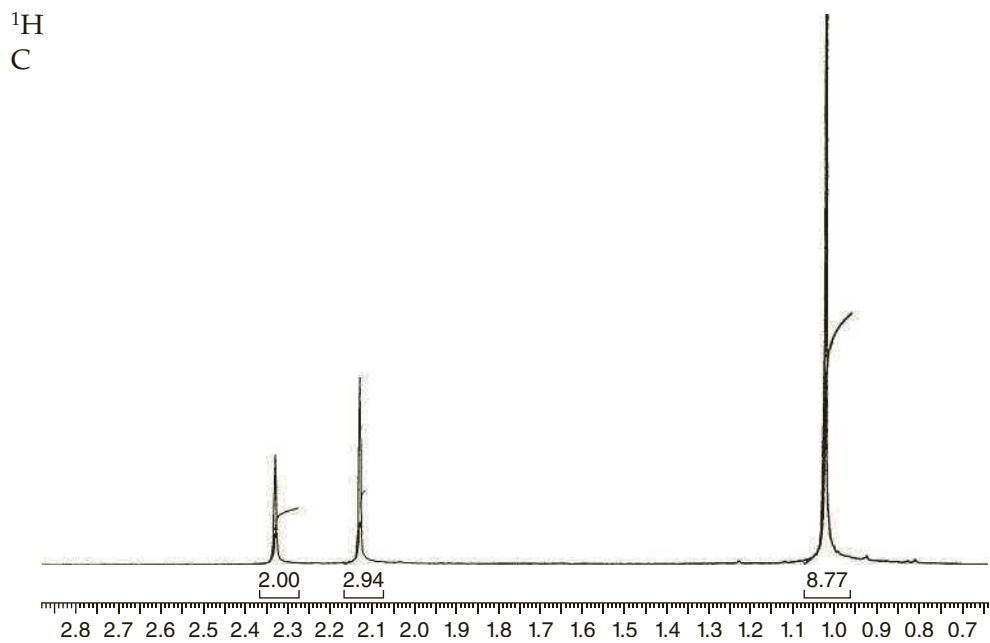

Mass Spectrometry

In its simplest form, the mass spectrometer performs three essential functions. First, molecules are bombarded by a stream of high-energy electrons, converting some of the molecules to positive ions. Because of their high energy, some of these ions **fragment**, or break apart into smaller ions. All of these ions are accelerated in an electric field. Second, the accelerated ions are separated according to their mass-to-charge ratio in a magnetic or electric field. Finally, the ions with a particular mass-to-charge ratio are detected by a device that is able to count the number of ions that strike it. The output of the detector is amplified and fed to a recorder. The trace from the recorder is a **mass spectrum**—a graph of the number of particles detected as a function of mass-to-charge ratio.

Ions are formed in an **ionization chamber**. The sample is introduced into the ionization chamber using a sample inlet system. In the ionization chamber, a heated **filament** emits a beam of high-energy electrons. The filament is heated to several thousand degrees Celsius. In normal operation, the electrons have an energy of about 70 electron-volts. These high-energy electrons strike a stream of molecules that has been admitted from the sample system and ionize the molecules in the sample stream by removing electrons from them. The molecules are thus converted into **radical-cations**.




The energy required to remove an electron from an atom or molecule is its **ionization potential**. The ionized molecules are accelerated and focused into a beam of rapidly moving ions by means of charged plates.

From the ionization chamber, the beam of ions passes through a short field-free region. From there, the beam enters the **mass analyzer**, where the ions are separated according to their mass-to-charge ratio.

The detector of most instruments consists of a counter that produces a current proportional to the number of ions that strike it. Electron multiplier circuits allow accurate measurement of the current from even a single ion striking the detector. The signal from the detector is fed to a **recorder**, which produces the actual mass spectrum.

28.1 The Mass Spectrum

The **mass spectrum** is a plot of ion abundance versus mass-to-charge (m/e) ratio. A typical mass spectrum is shown in Figure 28.1. The spectrum shown is that of dopamine, a substance that acts as a neurotransmitter in the central nervous system. The spectrum is displayed as a bar graph of percentage ion abundance (relative abundance) plotted against m/e .

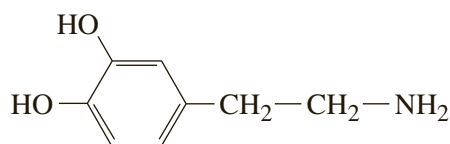

Dopamine

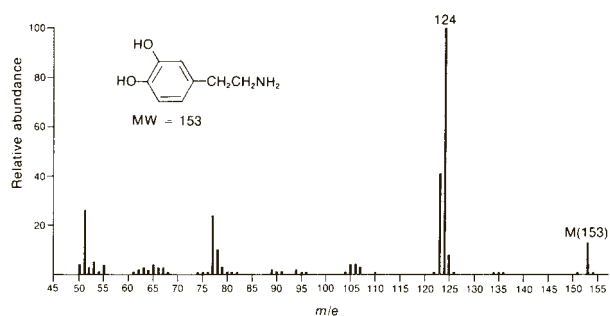


Figure 28.1
The mass spectrum of dopamine.

The most abundant ion formed in the ionization chamber gives rise to the tallest peak in the mass spectrum, called the **base peak**. For dopamine, the base peak appears at $m/e = 124$. The relative abundances of all the other peaks in the spectrum are reported as percentages of the abundance of the base peak.

The beam of electrons in the ionization chamber converts some of the sample molecules into positive ions. Removal of a single electron from a molecule yields an ion whose weight is the actual molecular weight of the original molecule. This ion is the **molecular ion**, frequently symbolized as M^+ .

The value of m/e at which the molecular ion appears on the mass spectrum, assuming that the ion has only one electron removed, gives the molecular weight of the original molecule. In the mass spectrum of dopamine, the molecular ion appears at $m/e = 153$, the molecular weight of dopamine. If you can identify the molecular ion peak in the mass spectrum, you can use the spectrum to determine the molecular weight of an unknown substance. If the presence of heavy isotopes is ignored for the moment, the molecular ion peak corresponds to the heaviest particle observed in the mass spectrum.

Molecules do not occur in nature as isotopically pure species. Virtually all atoms have heavier isotopes that occur in varying natural abundances. Hydrogen occurs largely as ^1H , but a few percent of hydrogen atoms occur as the isotope ^2H . Further, carbon normally occurs as ^{12}C , but a few percent of carbon atoms are the heavier isotope, ^{13}C . With the exception of fluorine, most other elements have a certain percentage of heavier isotopes that occur naturally. Peaks caused by ions bearing these heavier isotopes are also found in the mass spectrum. The relative abundances of these isotopic peaks are proportional to the abundances of the isotopes in nature. Most often, the isotopes occur at one or two mass units above the mass of the “normal” atom. Therefore, besides looking for the molecular ion (M^+) peak, you should also attempt to locate the $M+1$ and $M+2$ peaks. As will be demonstrated later, you can use the relative abundances of these $M+1$ and $M+2$ peaks to determine the molecular formula of the substance being studied.

The beam of electrons in the ionization chamber can produce the molecular ion. This beam also has sufficient energy to break some of the bonds in the molecule, producing a series of molecular fragments. Fragments that are positively charged are also accelerated in the ionization chamber, sent through the analyzer, detected, and recorded on the mass spectrum. These **fragment ion peaks** appear at m/e values corresponding to their individual masses. Very often, a fragment ion rather than the molecular ion will be the most abundant ion produced in the mass spectrum (the base peak).

A second means of producing fragment ions occurs with the molecular ion, which, once it is formed, is so unstable that it disintegrates before it can pass into the accelerating region of the ionization chamber. Lifetimes shorter than 10^{-5} seconds are typical in this type of fragmentation. Those fragments that are charged then appear as fragment ions in the mass spectrum. As a result of these fragmentation processes, the typical mass spectrum can be quite complex, containing many more peaks than the molecular ion and $M+1$ and $M+2$ peaks. Structural information about a substance can be determined by examining the fragmentation pattern in the mass spectrum. Fragmentation patterns are discussed further in Section 28.3.

28.2 Molecular Formula Determination

Mass spectrometry can be used to determine the molecular formulas of molecules that provide reasonably abundant molecular ions. Although there are at least two principal techniques for determining a molecular formula, only one will be described here.

The molecular formula of a substance can be determined through the use of **precise atomic masses**. High-resolution mass spectrometers are required for this method. Atoms are normally thought of as having integral atomic masses; for example, H = 1, C = 12, and O = 16. If you can determine atomic masses with sufficient precision, however, you find that the masses do not have values that are exactly integral. The mass of each atom actually differs from a whole mass number by a small fraction of a mass unit. The actual masses of some atoms are given in Table 28.1.

Depending on the atoms that are contained within a molecule, it is possible for particles of the same nominal mass to have slightly different measured masses when precise mass determinations are possible. To illustrate, a molecule whose molecular weight is 60 could be C_3H_8O , $C_2H_8N_2$, $C_2H_4O_2$, or CH_4N_2O . The species have the following precise masses:

C_3H_8O	60.05754
$C_2H_8N_2$	60.06884
$C_2H_4O_2$	60.02112
CH_4N_2O	60.03242

Table 28.1 Precise masses of some common elements

<i>Element</i>	<i>Atomic Weight</i>	<i>Nuclide</i>	<i>Precise Mass</i>
Hydrogen	1.00797	1H	1.00783
		2H	2.01410
Carbon	12.01115	^{12}C	12.0000
		^{13}C	13.00336
Nitrogen	14.0067	^{14}N	14.0031
		^{15}N	15.0001
Oxygen	15.9994	^{16}O	15.9949
		^{17}O	16.9991
		^{18}O	17.9992
Fluorine	18.9984	^{19}F	18.9984
Silicon	28.086	^{28}Si	27.9769
		^{29}Si	28.9765
		^{30}Si	29.9738
Phosphorus	30.974	^{31}P	30.9738
Sulfur	32.064	^{32}S	31.9721
		^{33}S	32.9715
		^{34}S	33.9679
Chlorine	35.453	^{35}Cl	34.9689
		^{37}Cl	36.9659
Bromine	79.909	^{79}Br	78.9183
		^{81}Br	80.9163
Iodine	126.904	^{127}I	126.9045

Observing a molecular ion with a mass of 60.058 would establish that the unknown molecule was C_3H_8O . Distinguishing among these possibilities is well within the capability of a modern high-resolution instrument.

In another method, these four compounds may also be distinguished by differences in the relative intensities of their M , $M+1$, and $M+2$ peaks. The predicted intensities are either calculated by formula or looked up in tables. Details of this method may be found in the References at the end of this chapter.

28.3 Detecting Halogens


When chlorine or bromine is present in a molecule, the isotope peak that is two mass units heavier than the molecular ion (the $M+2$ peak) becomes very significant. The heavy isotope of each of these elements is two mass units heavier than the lighter isotope. The natural abundance of ^{37}Cl is 32.5% that of ^{35}Cl ; the natural abundance of ^{81}Br is 98.0% that of ^{79}Br . When these elements are present, the $M+2$ peak becomes intense, and the pattern is characteristic of the particular halogen present. If a compound contains two chlorine or bromine atoms, a distinct $M+4$ peak should be observed, as well as an intense $M+2$ peak. In these cases, you should exercise caution in identifying the molecular ion peak in a mass spectrum, but the pattern of peaks is characteristic of the nature of the halogen substitution in the molecule. Table 28.2 gives the relative intensities of isotope peaks for various combinations of bromine and chlorine atoms. The patterns of molecular ion and isotopic peaks observed with halogen substitution are shown in Figure 28.2. Examples of these patterns can be seen in the mass spectra of chloroethane (Figure 28.3) and bromoethane (Figure 28.4).

28.4 Fragmentation Patterns


When the molecule has been bombarded by high-energy electrons in the ionization chamber of a mass spectrometer, besides losing one electron to form an ion, the molecule also absorbs some of the energy transferred in the collision between the molecule and the incident electrons. This extra energy puts the molecular ion in an excited vibrational state. The vibrationally excited molecular ion is often unstable and may lose some of this extra energy by breaking apart into fragments. If the lifetime of an individual molecular ion is longer than 10^{-5} seconds, a peak corresponding to the molecular ion will be observed in the mass spectrum. Those molecular ions with lifetimes shorter than 10^{-5} seconds will break apart into fragments before they are accelerated within the ionization chamber. In such cases, peaks

Table 28.2 Relative intensities of isotope peaks for various combinations of bromine and chlorine

Halogen	M	$M+2$	$M+4$	$M+6$
Br	100	97.7	—	—
Br ₂	100	195.0	95.4	—
Br ₃	100	293.0	286.0	93.4
Cl	100	32.6	—	—
Cl ₂	100	65.3	10.6	—
Cl ₃	100	97.8	31.9	3.47
BrCl	100	130.0	31.9	—
Br ₂ Cl	100	228.0	159.0	31.2
BrCl ₂	100	163.0	74.4	10.4

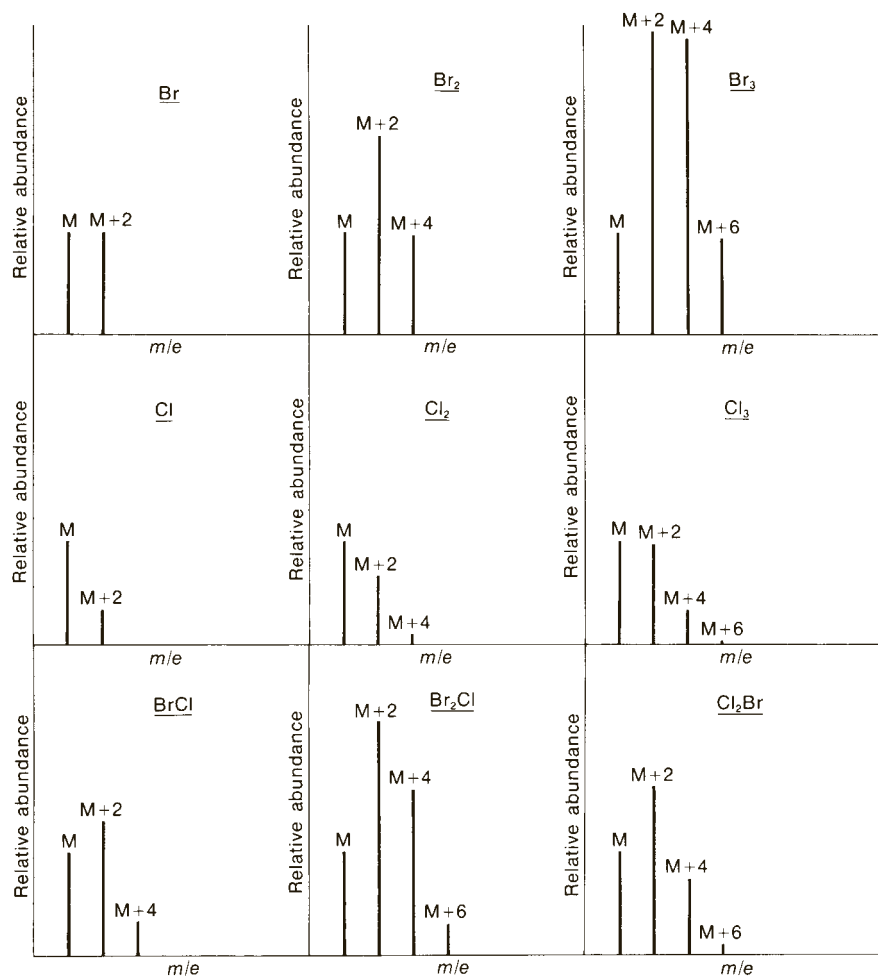**Figure 28.2**

Mass spectra expected for various combinations of bromine and chlorine.

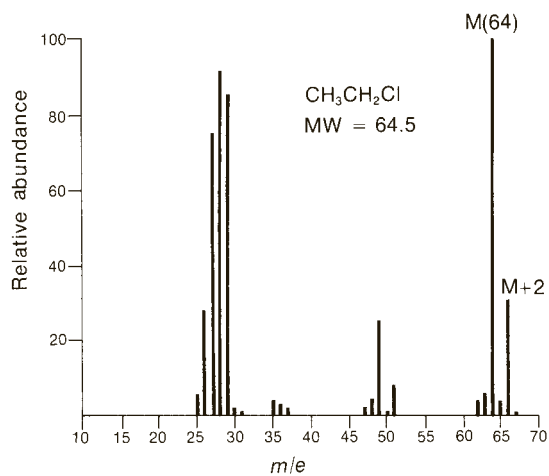**Figure 28.3**

The mass spectrum of chloroethane.

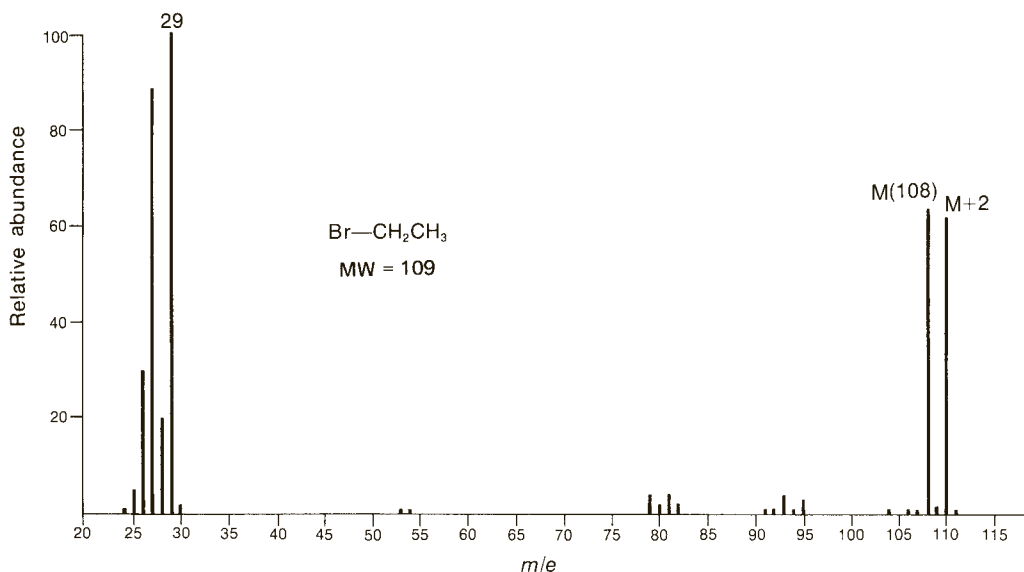


Figure 28.4
The mass spectrum of bromoethane.

corresponding to the mass-to-charge ratios for these fragments will also appear in the mass spectrum. For a given compound, not all the molecular ions formed by ionization have precisely the same lifetime. The ions have a range of lifetimes; some individual ions may have shorter lifetimes than others. As a result, peaks are usually observed arising from both the molecular ion and the fragment ions in a typical mass spectrum.

For most classes of compounds, the mode of fragmentation is somewhat characteristic. In many cases, it is possible to predict how a molecule will fragment. Remember that the ionization of the sample molecule forms a molecular ion that not only carries a positive charge but also has an unpaired electron. The molecular ion, then, is actually a **radical-cation**, and it contains an odd number of electrons. In the structural formulas that follow, the radical-cation is indicated by enclosing the structure in square brackets. The positive charge and the unshared electron are shown as superscripts.



When fragment ions form in the mass spectrometer, they almost always form by means of unimolecular processes. The pressure of the sample in the ionization chamber is too low to permit a significant number of bimolecular collisions. Those unimolecular processes that require the least energy will give rise to the most abundant fragment ions.


Fragment ions are cations. Much of the chemistry of these fragment ions can be explained in terms of what is known about carbocations in solution. For example, alkyl substitution stabilizes fragment ions (and promotes their formation) in much the same way that it stabilizes carbocations. Those fragmentation processes that lead to more stable ions will be favored over processes that lead to the formation of less stable ions.

Fragmentation often involves the loss of an electrically neutral fragment. The neutral fragment does not appear in the mass spectrum, but you can deduce its existence by noting the difference in masses of the fragment ion and the original molecular ion. Again, processes that lead to the formation of a more stable neutral fragment will be favored over those that lead to the formation of a less stable




neutral fragment. The loss of a stable neutral molecule, such as water, is commonly observed in the mass spectrometer.

A. Cleavage of One Bond

The most common mode of fragmentation involves the cleavage of one bond. In this process, the odd-electron molecular ion yields an odd-electron neutral fragment and an even-electron fragment ion. The neutral fragment that is lost is a **free radical**, whereas the ionic fragment is of the carbocation type. Cleavages that lead to the formation of more stable carbocations will be favored. Thus, the ease of fragmentation to form ions increases in the following order:

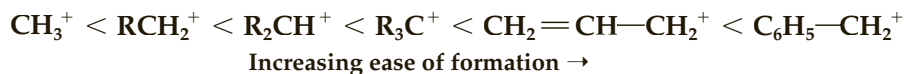



The following reactions show examples of fragmentation that take place with the cleavage of one bond:

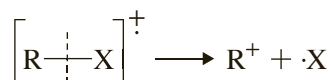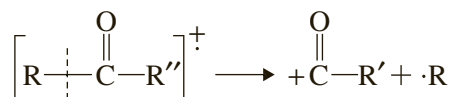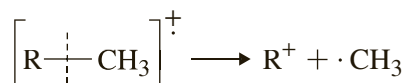

where X = halogen, OR,
SR, or NR₂, and where
R = H, alkyl, or aryl

B. Cleavage of Two Bonds

The next most important type of fragmentation involves the cleavage of two bonds. In this type of process, the odd-electron molecular ion yields an odd-electron fragment ion and an even-electron neutral fragment, usually a small, stable molecule. Examples of this type of cleavage are shown next:

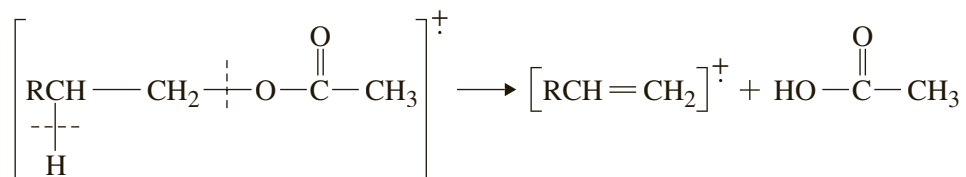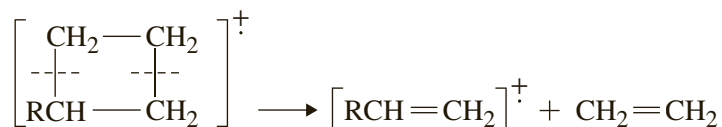


C. Other Cleavage Processes

In addition to the processes just mentioned, fragmentation reactions involving rearrangements, migrations of groups, and secondary fragmentations of fragment ions are also possible. These processes occur less often than the types of processes just described. Nevertheless, the pattern of molecular ion and fragment ion peaks observed in the typical mass spectrum is quite complex and unique for each particular molecule. As a result, the mass spectral pattern observed for a given substance can be compared with the mass spectra of known compounds as a means of identification. The mass spectrum is like a fingerprint. For a treatment of the specific modes of fragmentation characteristic of particular classes of compounds, refer to more advanced textbooks (see References at the end of this chapter). The unique appearance of the mass spectrum for a given compound is the basis for identifying the components of a mixture in the **gas chromatography–mass spectrometry (GC–MS)** technique (see Technique 22, Section 22.14). The mass spectrum of every component in a mixture is compared with standard spectra stored in the computer memory of the instrument. The printed output produced by a GC–MS instrument includes an identification based on the results of the computer matching of mass spectra.

28.5 Interpreted Mass Spectra

In this section, the mass spectra of some representative organic compounds are presented. The important fragment ion peaks in each mass spectrum are identified. In some of the examples, identification of the fragments is presented without explanation, although some interpretation is provided where an unusual or interesting process takes place. In the first example, that of butane, a more complete explanation of the symbolism used is offered.

Butane; C_4H_{10} , MW = 58 (Figure 28.5)

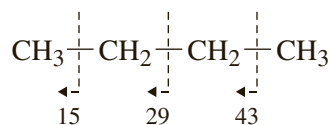

In the structural formula of butane, the dashed lines represent the location of bond-breaking processes that occur during fragmentation. In each case, the


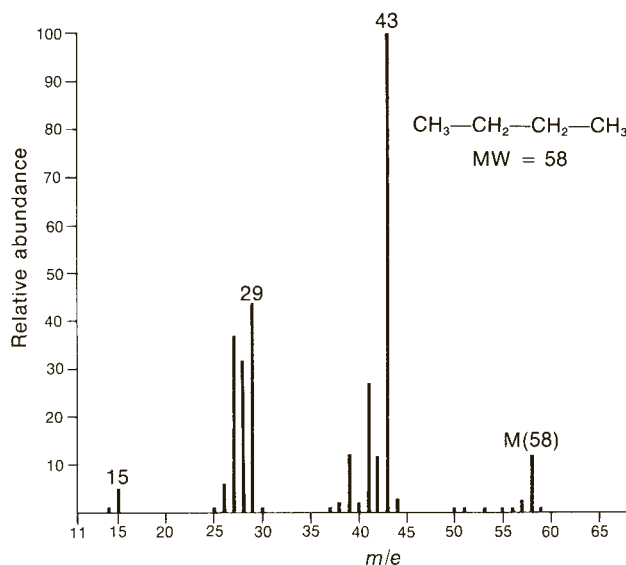



Figure 28.5
The mass spectrum of butane.

fragmentation process involves the breaking of one bond to yield a neutral radical and a cation. The arrows point toward the fragment that bears the positive charge. This positive fragment is the ion that appears in the mass spectrum. The mass of the fragment ion is indicated beneath the arrow.


The mass spectrum shows the molecular ion at $m/e = 58$. Breaking of the C1—C2 bond yields a three-carbon fragment with a mass of 43.

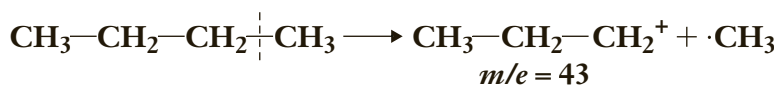

Cleavage of the central bond yields an ethyl cation, with a mass of 29.

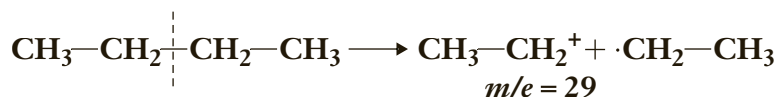

The terminal bond can also break to yield a methyl cation, which has a mass of 15.

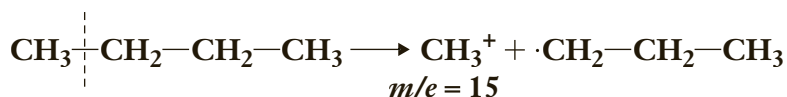

Each of these fragments appears in the mass spectrum of butane and has been identified.

2,2,4-Trimethylpentane; C_8H_{18} , MW = 114 (Figure 28.6)



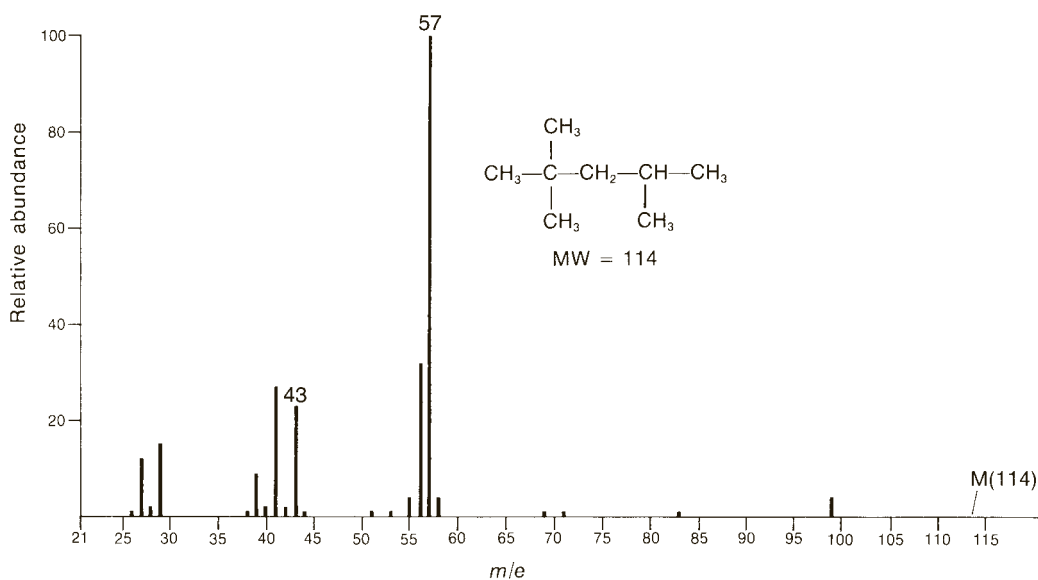


Figure 28.6
The mass spectrum of 2,2,4-trimethylpentane (“isooctane”).


Notice that in the case of 2,2,4-trimethylpentane, by far the most abundant fragment is the *tert*-butyl cation ($m/e = 57$). This result is not surprising when one considers that the *tert*-butyl cation is a particularly stable carbocation.

Cyclopentane; C_5H_{10} , MW = 70 (Figure 28.7)

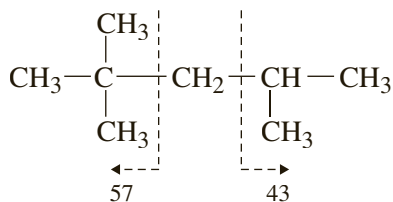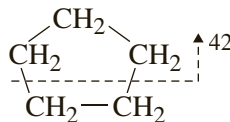


In the case of cyclopentane, the most abundant fragment results from the simultaneous cleavage of two bonds. This mode of fragmentation eliminates a neutral molecule of ethene (MW = 28) and results in the formation of a cation at $m/e = 42$.


1-Butene; C_4H_8 , MW = 56 (Figure 28.8)

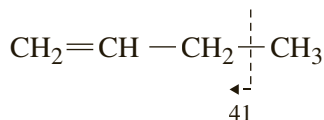

An important fragment in the mass spectra of alkenes is the allyl cation ($m/e = 41$). This cation is particularly stable due to resonance.

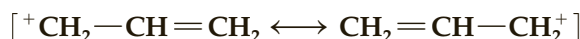

Toluene; C_7H_8 , MW = 92 (Figure 28.9)

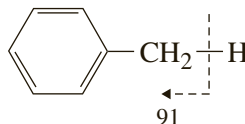

When an alkyl group is attached to a benzene ring, preferential fragmentation occurs at a benzylic position to form a fragment ion of the formula $C_7H_7^+$ ($m/e = 91$). In the mass spectrum of toluene, loss of hydrogen from the molecular ion gives a strong peak at $m/e = 91$. Although it may be expected that this fragment ion peak is due to the benzyl carbocation, evidence suggests the benzyl carbocation

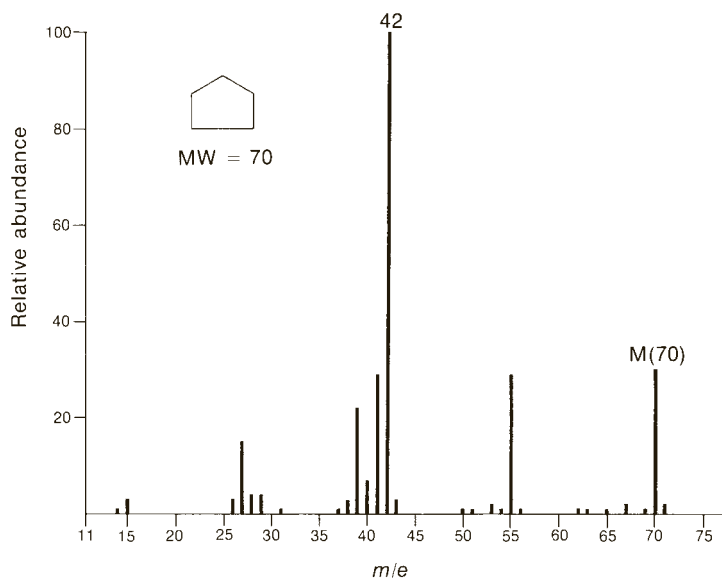


Figure 28.7
The mass spectrum of cyclopentane.

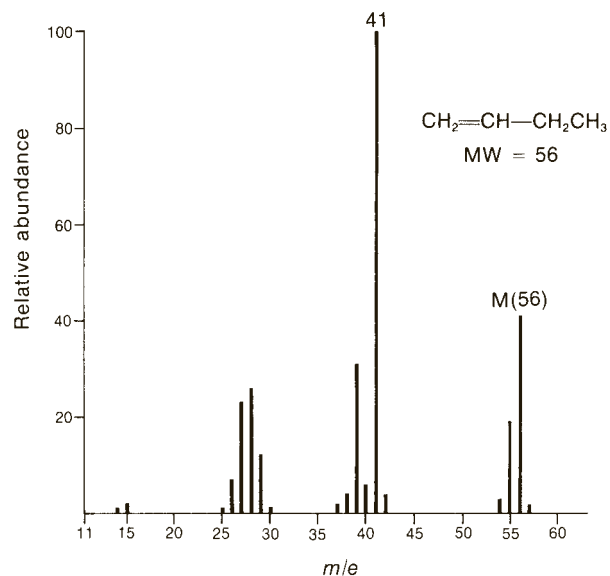


Figure 28.8
The mass spectrum of 1-butene.

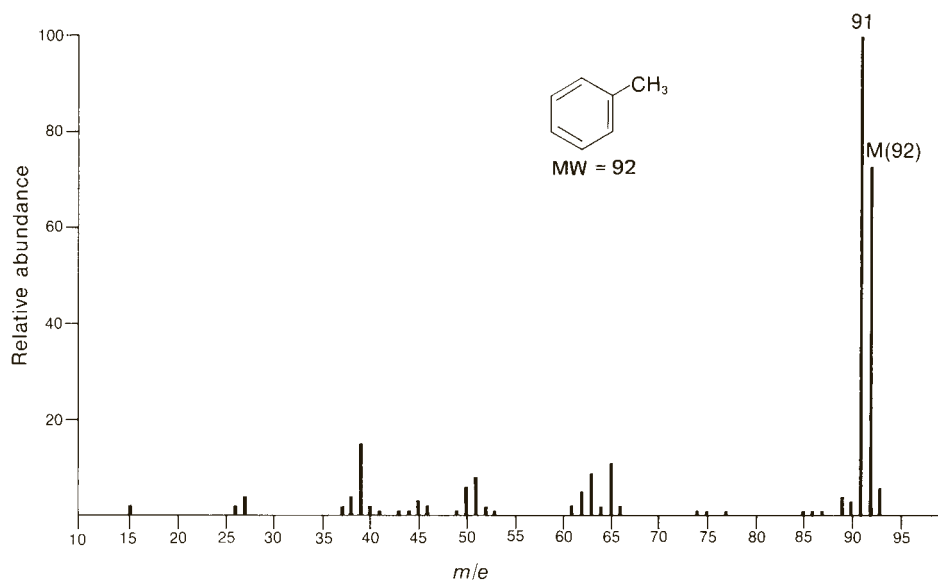
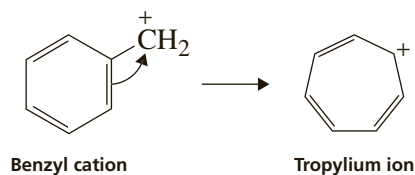


Figure 28.9
The mass spectrum of toluene.

actually rearranges to form the **tropylium ion**. Isotope-labeling experiments tend to confirm the formation of the tropylium ion. The tropylium ion is a seven-carbon ring system that contains six electrons in π -molecular orbitals and hence is resonance stabilized in a manner similar to that observed in benzene.





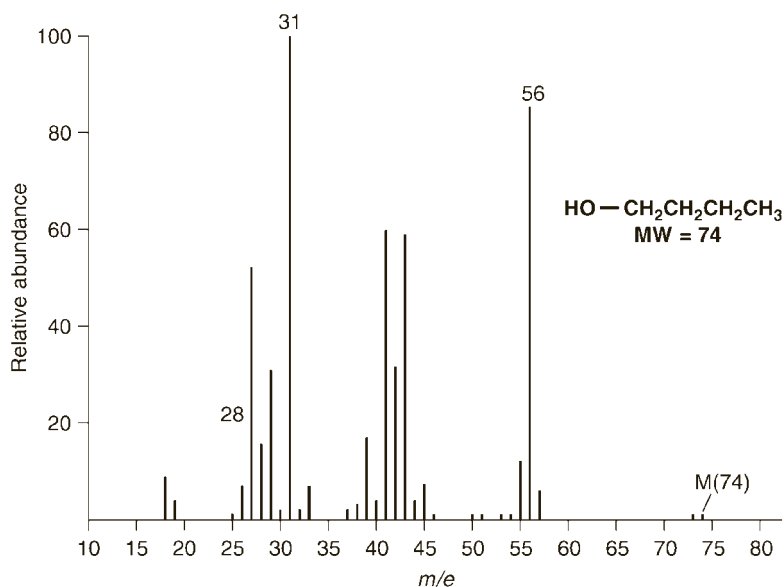


Figure 28.10
The mass spectrum of 1-butanol.


1-Butanol; $\text{C}_4\text{H}_{10}\text{O}$, MW = 74 (Figure 28.10)

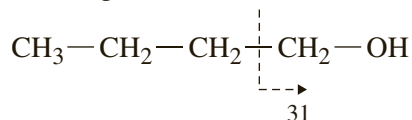


The most important fragmentation reaction for alcohols is loss of an alkyl group:

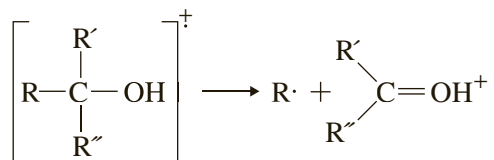


The largest alkyl group is the one that is lost most readily. In the spectrum of 1-butanol, the intense peak at $m/e = 31$ is due to the loss of a propyl group to form

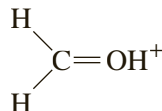

A second common mode of fragmentation involves dehydration. Loss of a molecule of water from 1-butanol leaves a cation of mass 56.

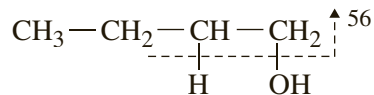

Benzaldehyde; $\text{C}_7\text{H}_6\text{O}$, MW = 106 (Figure 28.11)

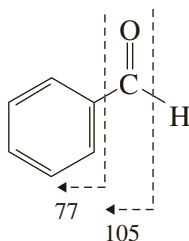


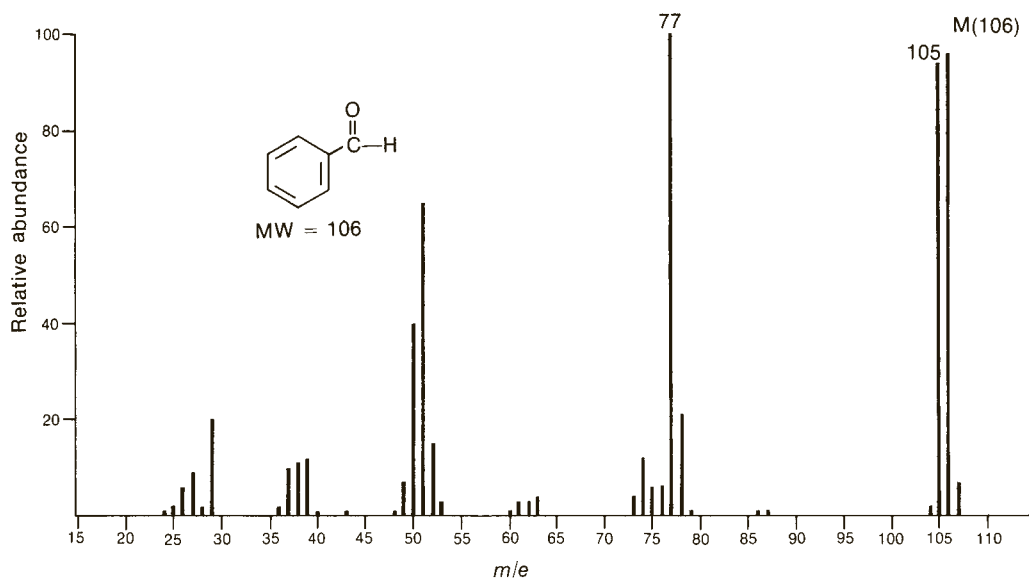



Figure 28.11
The mass spectrum of benzaldehyde.

The loss of a hydrogen atom from an aldehyde is a favorable process. The resulting fragment ion is a benzoyl cation, a particularly stable type of carbocation.

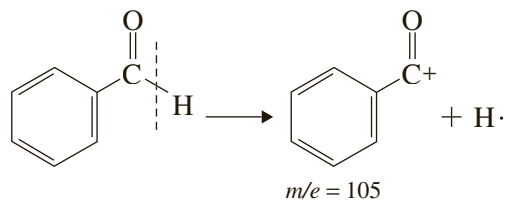


Loss of the entire aldehyde functional group leaves a phenyl cation. This ion can be seen in the spectrum of an m/e value of 77.

2-Butanone; C_4H_8O , MW = 72 (Figure 28.12)

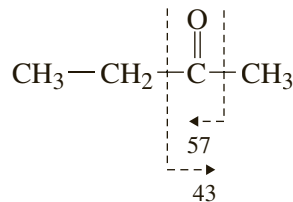

If the methyl group is lost as a neutral fragment, the resulting cation, an **acylium ion**, has an m/e value of 57. If the ethyl group is lost, the resulting acylium ion appears at an m/e value of 43.

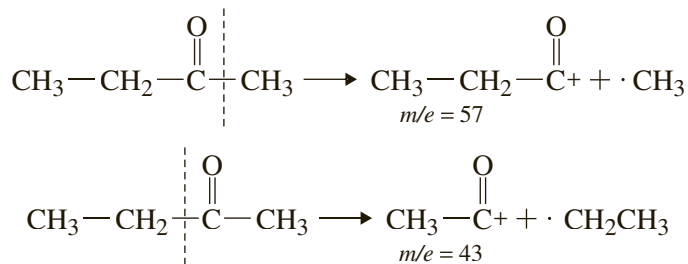

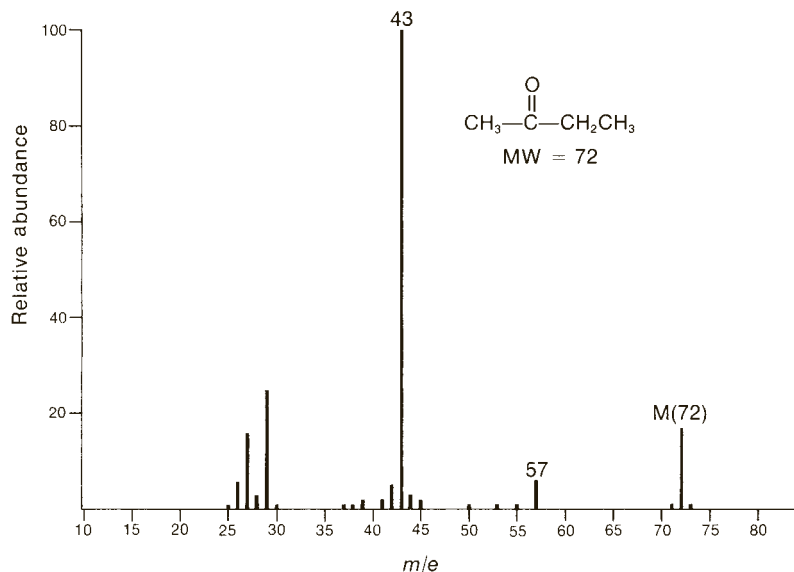


Figure 28.12
The mass spectrum of 2-butanone.



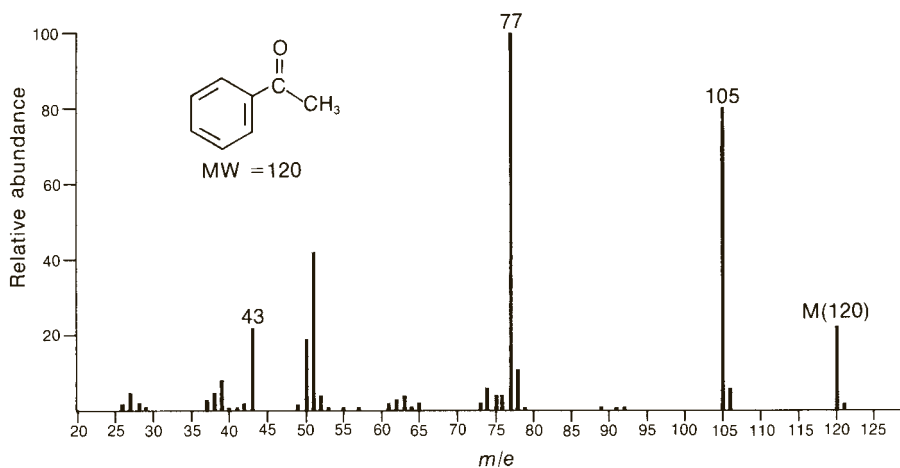


Figure 28.13
The mass spectrum of acetophenone.


Acetophenone; $\text{C}_8\text{H}_8\text{O}$, MW = 120 (Figure 28.13)

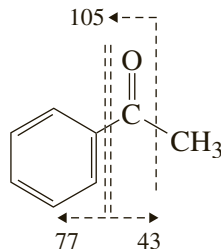

Aromatic ketones undergo α -cleavage to lose the alkyl group and form the benzoyl cation ($m/e = 105$). This ion subsequently loses carbon monoxide to form the phenyl cation ($m/e = 77$). Aromatic ketones also undergo α -cleavage on the other side of the carbonyl group, forming an alkyl acylium ion. In the case of acetophenone, this ion appears at an m/e value of 43.

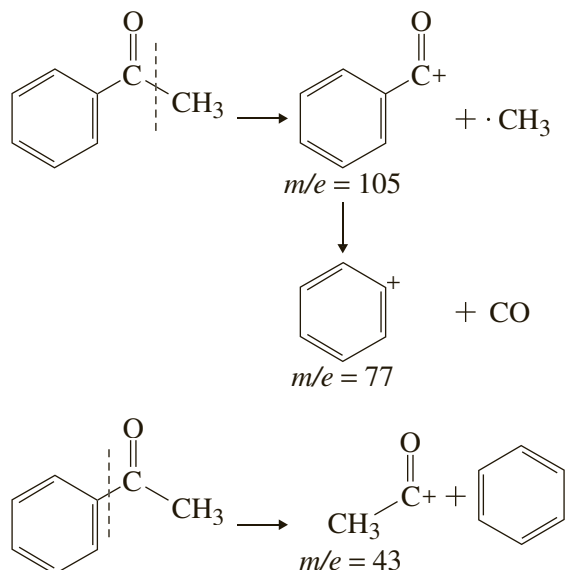

Propanoic acid; C₃H₆O₂, MW = 74 (Figure 28.14)

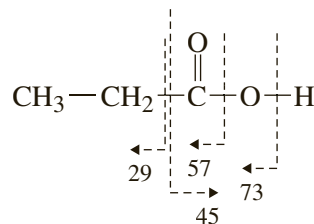

With short-chain carboxylic acids, the loss of OH and COOH through α -cleavage on either side of the C = O group may be observed. In the mass spectrum of propanoic acid, loss of OH gives rise to a peak at $m/e = 57$. Loss of COOH gives rise to a peak at $m/e = 29$. Loss of the alkyl group as a free radical, leaving the COOH⁺ ion ($m/e = 45$), also occurs. The intense peak at $m/e = 28$ is due to additional fragmentation of the ethyl portion of the acid molecule.

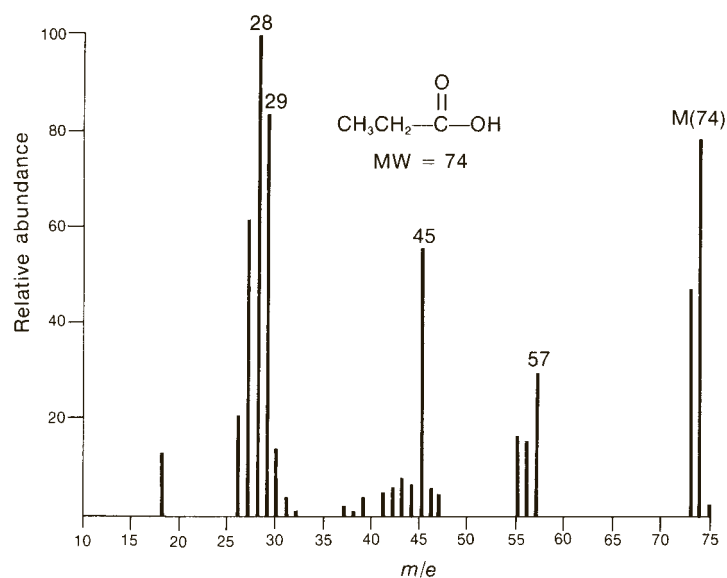


Figure 28.14
The mass spectrum of propanoic acid.


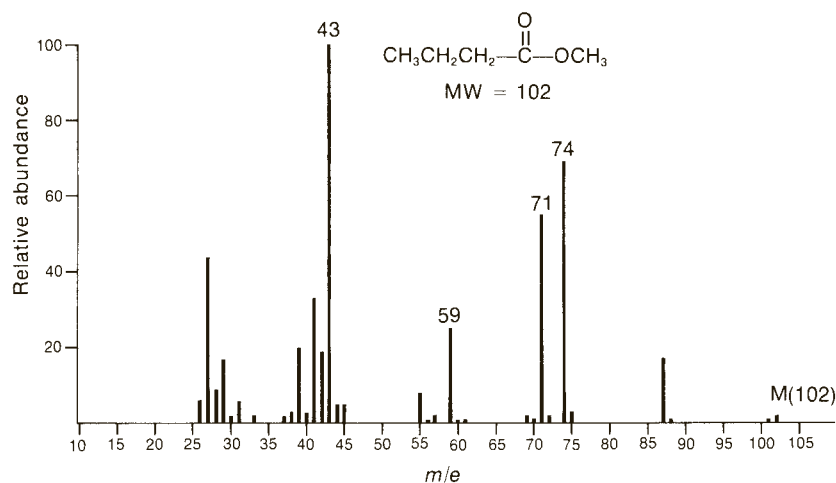



Figure 28.15
The mass spectrum of methyl butanoate.

Methyl butanoate; $\text{C}_5\text{H}_{10}\text{O}_2$, MW = 102 (Figure 28.15)

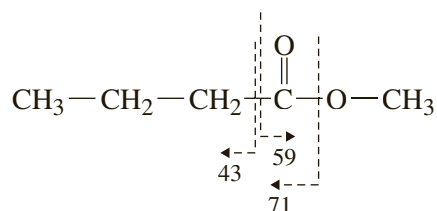

The most important of the α -cleavage reactions involves the loss of the alkoxy group from the ester to form the corresponding acylium ion, RCO^+ . The acylium ion peak appears at $m/e = 71$ in the mass spectrum of methyl butanoate. A second important peak results from the loss of the alkyl group from the acyl portion of the ester molecule, leaving a fragment $\text{CH}_3-\text{O}-\text{C}=\text{O}^+$ that appears at $m/e = 59$. Loss of the carboxylate function group to leave the alkyl group as a cation gives rise to a peak at $m/e = 43$. The intense peak at $m/e = 74$ results from a rearrangement process (see Section 28.6).

1-Bromohexane; $\text{C}_6\text{H}_{13}\text{Br}$, MW = 165 (Figure 28.16)

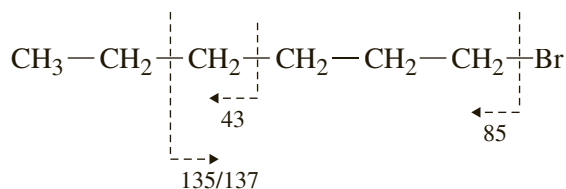

The most interesting characteristic of the mass spectrum of 1-bromohexane is the presence of the doublet in the molecular ion. These two peaks, of equal height, separated by two mass units, are strong evidence that bromine is present in the substance. Notice also that loss of the terminal ethyl group yields a fragment ion that still contains bromine ($m/e = 135$ and 137). The presence of the doublet demonstrates that this fragment contains bromine.

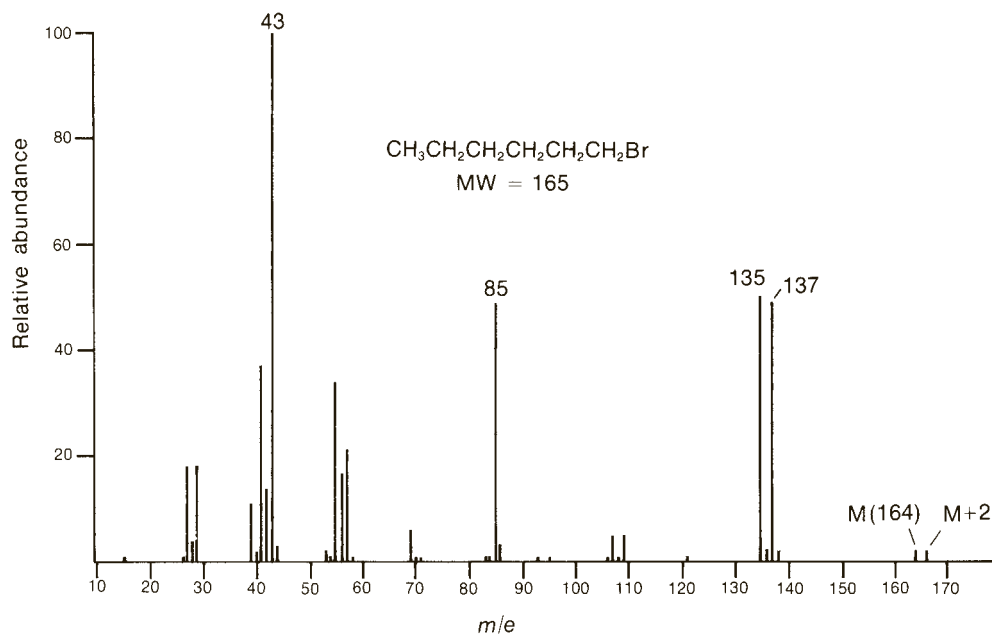



Figure 28.16
The mass spectrum of 1-bromohexane.

28.6 Rearrangement Reactions

Because the fragment ions that are detected in a mass spectrum are cations, we can expect that these ions will exhibit behavior we are accustomed to associate with carbocations. It is well known that carbocations are prone to rearrangement reactions, converting a less stable carbocation into a more stable one. These types of rearrangements are also observed in the mass spectrum. If the abundance of a cation is especially high, it is assumed that a rearrangement to yield a longer-lived cation must have occurred.

Other types of rearrangements are also known. An example of a rearrangement that is not normally observed in solution chemistry is the rearrangement of a benzyl cation to a tropylium ion. This rearrangement is seen in the mass spectrum of toluene (Figure 28.9).

A particular type of rearrangement process that is unique to mass spectrometry is the **McLafferty rearrangement**. This type of rearrangement occurs when an alkyl chain of at least three carbons in length is attached to an energy-absorbing structure such as a phenyl or carbonyl group that can accept the transfer of a hydrogen ion. The mass spectrum of methyl butanoate (Figure 28.15) contains a prominent peak at $m/e = 74$. This peak arises from a McLafferty rearrangement of the molecular ion.

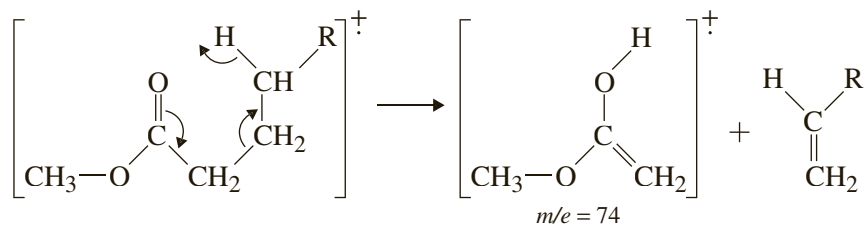

REFERENCES

- Beynon, J. H. *Mass Spectrometry and Its Applications to Organic Chemistry*. Elsevier: Amsterdam, **1960**.
- Biemann, K. *Mass Spectrometry: Organic Chemical Applications*. McGraw-Hill: New York, **1962**.
- Budzikiewicz, H.; Djerassi, C.; Williams, D. H. *Mass Spectrometry of Organic Compounds*. Holden-Day: San Francisco, **1967**.
- McLafferty, F. W.; Tureček, F. *Interpretation of Mass Spectra*, 4th ed. University Science Books: Mill Valley, CA, **1993**.
- Pavia, D.L.; Lampman, G. M.; Kriz, G.S.; and Vyvyan, J.R. *Introduction to Spectroscopy*, 5th ed., Cengage Learning: Stamford, CT, **2015**.
- Silverstein, R.M.; Webster, F.X.; and Kiemle, D. *Spectrometric Identification of Organic Compounds*, 7th ed., John Wiley: New York, **2005**.

Guide to the Chemical Literature

Often, you may need to go beyond the information contained in the typical organic chemistry textbook and use reference material in the library. At first glance, using library materials may seem formidable because of the numerous sources the library contains. If, however, you adopt a systematic approach, the task can prove rather useful. This description of various popular sources and an outline of logical steps to follow in the typical literature search should be helpful.

29.1 Locating Physical Constants: Handbooks

To find information on routine physical constants, such as melting points, boiling points, indices of refraction, and densities, you should first consider a handbook. Examples of suitable handbooks are

Aldrich Handbook of Fine Chemicals. Sigma-Aldrich: Milwaukee, WI, 2014–2016.

Haynes, W. M., ed. *CRC Handbook of Chemistry and Physics*, 94th ed., CRC Press: Boca Raton, FL, 2013.

O'Neil, M. S., ed. *The Merck Index*, 15th ed., Royal Society of Chemistry, London, 2013.

Speight, J., ed. *Lange's Handbook of Chemistry*, 70th Anniversary (16th) Edition. McGraw-Hill: New York, 2004.

Each of these references is discussed in detail in Technique 4. The *CRC Handbook* is the reference consulted most often because the book is so widely available. This handbook has been available, with periodic updates, for more than 100 years. It was first published in 1913. For those colleges, universities, and corporations that subscribe, an online version is now available. With the correct applet installed, the handbook can be searched using a smartphone. There is even a smartphone applet that will allow you to sketch the desired search structure. There are, however, distinct advantages to using the other handbooks. The *CRC Handbook* uses the *Chemical Abstracts* system of nomenclature in a format that requires you to identify the parent name; 3-methyl-1-butanol is listed as 1-butanol, 3-methyl. Other handbooks use simpler systems and may even include common names.

The *Merck Index*, which is also available to subscribers online, lists fewer compounds, but there is far more information provided for the ones listed. If the compound is a medicinal or natural product, this is the reference of choice. This handbook contains literature references for the isolation and synthesis of a compound, along with certain properties of medicinal interest, such as toxicity. *Lange's Handbook* and the *Aldrich Handbook* list compounds in alphabetical order; 3-methyl-1-butanol is listed as 3-methyl-1-butanol. The information in *Lange's Handbook* is available online to those who subscribe. Information about compounds in the *Aldrich Handbook* (including MSDS sheets) can be obtained from their website www.Sigma-Aldrich.com.

A more complete handbook that is usually housed in the library is

Buckingham, J., ed. *Dictionary of Organic Compounds*, 6th ed., Chapman & Hall/CRC Press: Boca Raton, FL, 1995.

This is a revised version of an earlier four-volume handbook edited by I. M. Heilbron and H. M. Bunbury. In its present form, it consists of nine volumes with two supplements.

29.2 General Synthetic Methods

Many standard introductory textbooks in organic chemistry provide tables that summarize most of the common reactions, including side reactions, for a given class of compounds. These books also describe alternative methods of preparing compounds.

Brown, W. H.; Foote, C. S.; Iverson, B. L.; Anslyn, E. *Organic Chemistry*, 7th ed., Brooks/Cole: Pacific Grove, CA, 2013.

Bruice, P. Y. *Organic Chemistry*, 7th ed., Prentice Hall: New York, 2013.

Carey, F. A.; Giuliano, R. M. *Organic Chemistry*, 8th ed., McGraw-Hill: New York, 2011.

Ege, S. *Organic Chemistry*, 5th ed., Houghton-Mifflin: Boston, 2004.

Fessenden, R. J.; Fessenden, J. S. *Organic Chemistry*, 6th ed., Brooks/Cole: Pacific Grove, CA, 1998.

Fox, M. A.; Whitesell, J. K. *Organic Chemistry*, 3rd ed., Jones & Bartlett: Boston, 2004.

Hornback, J. *Organic Chemistry*. 2nd ed., Brooks/Cole: Pacific Grove, CA, 2006.

Jones, M., Jr.; Fleming, S. A. *Organic Chemistry*, 5th ed., W. W. Norton: New York, 2014.

Loudon, G. M. *Organic Chemistry*, 8th ed., Benjamin/Cummings: Menlo Park, CA, 2012.

McMurry, J. *Organic Chemistry*, 8th ed., Brooks/Cole: Pacific Grove, CA, 2011.

Morrison, R. T.; Boyd, R. N. *Organic Chemistry*, 8th ed., Prentice Hall: Englewood Cliffs, NJ, 2012.

Smith, M. B.; March, J. *Advanced Organic Chemistry*, 7th ed., John Wiley & Sons: New York, 2013.

Solomons, T. W. G.; Fryhle, C.; Snyder, S. A. *Organic Chemistry*, 11th ed., John Wiley & Sons: New York, 2013.

Streitwieser, A.; Heathcock, C. H.; Kosower, E. M. *Introduction to Organic Chemistry*, Rev. 4th ed., Prentice Hall: New York, 1998.

Vollhardt, K. P. C.; Schore, N. E. *Organic Chemistry*, 6th ed., W. H. Freeman: New York, 2010.

Wade, L. G., Jr. *Organic Chemistry*, 8th ed., Prentice Hall: Englewood Cliffs, NJ, 2012.

29.3 Searching the Chemical Literature

If the information you are seeking is not available in any of the handbooks mentioned in Section 29.1 or if you are searching for more detailed information than they can provide, then a proper literature search is in order. Although an examination of standard textbooks can provide some help, you often must use all the resources of the library, including journals, reference collections, and abstracts. The following sections outline how the various types of sources should be used and what sort of information can be obtained from them.

The methods discussed for searching the literature use mainly printed materials. Modern search methods also make use of computerized databases and are discussed in Section 29.11. These are vast collections of data and bibliographic materials that can be scanned rapidly from remote computer terminals. Although computerized searching is widely available, it may not be readily accessible to

undergraduate students. The following references provide excellent introductions to the literature of organic chemistry:

- Carr, C. Teaching and Using Chemical Information. *Journal of Chemical Education*, **1993**, 719.
- Maizell, R. E. *How to Find Chemical Information*, 3rd ed., John Wiley & Sons: New York, **1998**.
- Smith, M. B.; March, J. *Advanced Organic Chemistry*, 7th ed., John Wiley & Sons: New York, **2013**.
- Somerville, A. N. Information Sources for Organic Chemistry, 1: Searching by Name Reaction and Reaction Type. *Journal of Chemical Education*, **1991**, 553.
- Somerville, A. N. Information Sources for Organic Chemistry, 2: Searching by Functional Group. *Journal of Chemical Education*, **1991**, 842.
- Somerville, A. N. Information Sources for Organic Chemistry, 3: Searching by Reagent. *Journal of Chemical Education*, **1992**, 379.
- Wiggins, G. *Chemical Information Sources*. McGraw-Hill: New York, **1991**. Integrates printed materials and computer sources of information.

29.4 Collections of Spectra

Collections of infrared, nuclear magnetic resonance, and mass spectra can be found in the following catalog of spectra:

- Cornu, A.; Massot, R. *Compilation of Mass Spectral Data*, 2nd ed., Heyden and Sons: London, **1975**.
- High-Resolution NMR Spectra Catalog*. Varian Associates: Palo Alto, CA, Vol. 1, **1962**; Vol. 2, **1963**.
- Johnson, L. F.; Jankowski, W. C. *Carbon-13 NMR Spectra*. John Wiley & Sons: New York, **1972**.
- Pouchert, C. J. *Aldrich Library of Infrared Spectra*, 3rd ed., Aldrich Chemical Co.: Milwaukee, WI, **1981**.
- Pouchert, C. J. *Aldrich Library of FT-IR Spectra*, 2nd ed., Aldrich Chemical Co.: Milwaukee, WI, **1997**.
- Pouchert, C. J. *Aldrich Library of NMR Spectra*, 2nd ed., Aldrich Chemical Co.: Milwaukee, WI, **1983**.
- Pouchert, C. J.; Behnke, J. *Aldrich Library of ¹³C and ¹H FT NMR Spectra*. Aldrich Chemical Co.: Milwaukee, WI, **1993**.
- Sadtler Standard Spectra*. Sadtler Research Laboratories: Philadelphia Continuing collection.
- Stenhagen, E.; Abrahamsson, S.; McLafferty, F. W. *Registry of Mass Spectral Data*, 4 vols. Wiley-Interscience: New York, **1974**.

The American Petroleum Institute has also published collections of infrared, nuclear magnetic resonance, and mass spectra.

29.5 Advanced Textbooks

Much information about synthetic methods, reaction mechanisms, and reactions of organic compounds is available in any of the many advanced textbooks in organic chemistry. Examples of such books are

- Carey, F. A.; Sundberg, R. J. *Advanced Organic Chemistry. Part A. Structure and Mechanisms; Part B. Reactions and Synthesis*, 5th ed. Kluwer Academic: New York, **2008**.

- Carruthers, W. *Some Modern Methods of Organic Synthesis*, 4th ed., Cambridge University Press: Cambridge, UK, 2004.
- Corey, E. J.; Cheng, X-M. *The Logic of Chemical Synthesis*. John Wiley & Sons: New York, 1995.
- Fieser, L. F.; Fieser, M. *Advanced Organic Chemistry*. Reinhold: New York, 1961.
- Finar, I. L. *Organic Chemistry*, 6th ed., Longman Group: London, 1986.
- House, H. O. *Modern Synthetic Reactions*, 2nd ed., W. H. Benjamin: Menlo Park, CA, 1972.
- Noller, C. R. *Chemistry of Organic Compounds*, 3rd ed., W. B. Saunders: Philadelphia, 1965.
- Smith, M. B. *Organic Synthesis*, 2nd ed., McGraw-Hill: New York, 2002.
- Smith, M. B.; March, J. *Advanced Organic Chemistry*, 7th ed., John Wiley & Sons: New York, 2013.
- Stowell, J. C. *Intermediate Organic Chemistry*, 2nd ed., John Wiley & Sons: New York, 1993.
- Warren, S.; Wyatt, P. *Organic Synthesis: The Disconnection Approach*. 2nd ed., John Wiley & Sons: New York, 2009.
- Zweifel, G. S.; Nantz, M. H. *Modern Organic Synthesis*. W. H. Freeman and Company: New York, 2007.

These books often contain references to original papers in the literature for students wanting to follow the subject further. Consequently, you obtain not only a review of the subject from such a textbook but also a key reference that is helpful toward a more extensive literature search. The textbook by Smith and March is particularly useful for this purpose.

29.6 Specific Synthetic Methods

Anyone interested in locating information about a particular method of synthesizing a compound should first consult one of the many general textbooks on the subject. Useful ones are:

- Anand, N.; Bindra, J. S.; Ranganathan, S. *Art in Organic Synthesis*, 2nd ed., John Wiley & Sons: New York, 1988.
- Barton, D.; Ollis, W. D., eds. *Comprehensive Organic Chemistry*, 6 vols. Pergamon Press: Oxford, 1979.
- Buehler, C. A.; Pearson, D. E. *Survey of Organic Synthesis*. Wiley-Interscience: New York, 1970, 2 vols., 1977.
- Carey, F. A.; Sundberg, R. J. *Advanced Organic Chemistry. Part B. Reactions and Synthesis*, 5th ed., Kluwer: New York, 2008.
- Compendium of Organic Synthetic Methods*. Wiley-Interscience: New York, 1971–2009. This is a continuing series, now in 12 volumes.
- Fieser, L. F.; Fieser, M. *Reagents for Organic Synthesis*. Wiley-Interscience: New York, 1967–2013. This is a continuing series, now in 27 volumes.
- Greene, T. W.; Wuts, P. G. M. *Protective Groups in Organic Synthesis*, 4th ed., John Wiley & Sons: New York, 2007.
- House, H. O. *Modern Synthetic Reactions*, 2nd ed., W. H. Benjamin: Menlo Park, CA, 1972.
- Larock, R. C. *Comprehensive Organic Transformations*, 2nd ed., Wiley-VCH: New York, 1999.

Mundy, B. P.; Eller, M. G. *Name Reactions and Reagents in Organic Synthesis*. 2nd ed., John Wiley & Sons: New York, **2005**.

Patai, S., ed. *The Chemistry of the Functional Groups*. Interscience, 1964–present: London, **2005**. This series consists of many volumes, each one specializing in a particular functional group.

Smith, M. B.; March, J. *Advanced Organic Chemistry*, 7th ed., John Wiley & Sons: New York, **2013**.

Trost, B. M.; Fleming, I. *Comprehensive Organic Synthesis*. Pergamon/Elsevier Science: Amsterdam, **1992**. This series consists of 9 volumes plus supplements.

Vogel, A. I. *Vogel's Textbook of Practical Organic Chemistry, Including Qualitative Organic Analysis*, 5th ed., Longman Group: London, **1989**. Revised by members of the School of Chemistry, Thames Polytechnic.

Wagner, R. B.; Zook, H. D. *Synthetic Organic Chemistry*. John Wiley & Sons: New York, **1956**.

Wang, Z. *Comprehensive Organic Name Reactions and Reagents*. John Wiley: New York, **2009**.

More specific information, including actual reaction conditions, exists in collections specializing in organic synthetic methods. The most important of these are:

Organic Syntheses. John Wiley & Sons: New York, **1921**–present. Published annually.

Organic Syntheses, Collective Volumes. John Wiley & Sons: New York, **1941–2014**.

Vol. 1, 1941, Annual Volumes 1–9	Vol. 7, 1990, Annual Volumes 60–64
Vol. 2, 1943, Annual Volumes 10–19	Vol. 8, 1993, Annual Volumes 65–69
Vol. 3, 1955, Annual Volumes 20–29	Vol. 9, 1998, Annual Volumes 70–74
Vol. 4, 1963, Annual Volumes 30–39	Vol. 10, 2004, Annual Volumes 75–79
Vol. 5, 1973, Annual Volumes 40–49	Vol. 11, 2010, Annual Volumes 80–84
Vol. 6, 1988, Annual Volumes 50–59	Vol. 12, 2014, Annual Volumes 85–89

It is much more convenient to use the collective volumes where the earlier annual volumes of *Organic Syntheses* are combined in groups of 9 or 10 in the first six collective volumes (Volumes 1–6), and then in groups of 5 for the next four volumes (Volumes 7–10). Useful indices are included at the end of each of the collective volumes that classify methods according to the type of reaction, type of compound prepared, formula of compound prepared, preparation or purification of solvents and reagents, and use of various types of specialized apparatus.

The main advantage of using one of the *Organic Syntheses* procedures is that they have been tested to make sure that they work as written. Often, an organic chemist will adapt one of these tested procedures to the preparation of another compound. One of the features of the advanced organic textbook by Smith and March is that it includes references to specific preparative methods contained in *Organic Syntheses*.

More advanced material on organic chemical reactions and synthetic methods may be found in any one of a number of annual publications that review the original literature and summarize it. Examples include

Advances in Organic Chemistry: Methods and Results. John Wiley & Sons: New York, **1960**–present.

Annual Reports in Organic Synthesis. Academic Press: Orlando, FL, **1985**–present.

Annual Reports of the Chemical Society, Section B. Chemical Society: London, **1905**–present. Specifically, the section Synthetic Methods.

Organic Reactions. John Wiley & Sons: New York, **1942**–present.

Progress in Organic Chemistry. John Wiley & Sons: New York, 1952–1973.

Each of these publications contains a great many citations to the appropriate articles in the original literature.

29.7 Advanced Laboratory Techniques

The student who is interested in reading about techniques more advanced than those described in this textbook, or in more complete descriptions of techniques, should consult one of the advanced textbooks specializing in organic laboratory techniques. Besides focusing on apparatus construction and the performance of complex reactions, these books provide advice on purifying reagents and solvents. Useful sources of information on organic laboratory techniques include:

Bates, R. B.; Schaefer, J. P. *Research Techniques in Organic Chemistry*. Prentice Hall: Englewood Cliffs, NJ, 1971.

Krubsack, A. J. *Experimental Organic Chemistry*. Allyn & Bacon: Boston, 1973.

Leonard, J.; Lygo, B.; Procter, G. *Advanced Practical Organic Chemistry*, 2nd ed., Chapman & Hall: London, 1995.

Monson, R. S. *Advanced Organic Synthesis: Methods and Techniques*. Academic Press: New York, 1971.

Pirrung, M. C. *The Synthetic Organic Chemist's Companion*. John Wiley: New York, 2009.

Techniques of Chemistry. John Wiley & Sons: New York, 1970–present. Currently over 25 volumes. The successor to *Technique of Organic Chemistry*, this series covers experimental methods of chemistry, such as purification of solvents, spectral methods, and kinetic methods.

Weissberger, A., et al., eds. *Technique of Organic Chemistry*, 3rd ed., 14 vols. Wiley-Interscience: New York, 1959–1969.

Wiberg, K. B. *Laboratory Technique in Organic Chemistry*. McGraw-Hill: New York, 1960.

Numerous works and some general textbooks specialize in particular techniques. The preceding list is representative only of the most common books in this category. The following books deal specifically with microscale and semimicroscale techniques.

Cheronis, N. D.; Micro and Semimicro Methods. In A. Weissberger, ed., *Technique of Organic Chemistry*, Vol. 6. Wiley-Interscience: New York, 1954.

Cheronis, N. D.; Ma, T. S. *Organic Functional Group Analysis by Micro and Semimicro Methods*. Wiley-Interscience: New York, 1964.

Ma, T. S.; Horak, V. *Microscale Manipulations in Chemistry*. Wiley-Interscience: New York, 1976.

29.8 Reaction Mechanisms

As with the case of locating information on synthetic methods, you can obtain a great deal of information about reaction mechanisms by consulting one of the common textbooks on physical organic chemistry. The textbooks listed here provide a general description of mechanisms, but they do not contain specific literature citations. Very general textbooks include:

Bruckner, R. *Advanced Organic Chemistry: Reaction Mechanisms*. Academic Press: New York, 2001.

Miller, A.; Solomon, P. *Writing Reaction Mechanisms in Organic Chemistry*, 2nd ed., Academic Press: San Diego, CA, 1999.

Sykes, P. *A Primer to Mechanisms in Organic Chemistry*. Benjamin/Cummings: Menlo Park, CA, 1995.

More advanced textbooks include

Carey, F. A.; Sundberg, R. J. *Advanced Organic Chemistry. Part A. Structure and Mechanisms*, 5th ed., Kluwer: New York, 2008.

Hammett, L. P. *Physical Organic Chemistry: Reaction Rates, Equilibria, and Mechanisms*, 2nd ed., McGraw-Hill: New York, 1970.

Hine, J. *Physical Organic Chemistry*, 2nd ed., McGraw-Hill: New York, 1962.

Ingold, C. K. *Structure and Mechanism in Organic Chemistry*, 2nd ed., Cornell University Press: Ithaca, NY, 1969.

Isaacs, N. S. *Physical Organic Chemistry*, 2nd ed., John Wiley & Sons: New York, 1995.

Jones, R. A. Y. *Physical and Mechanistic Organic Chemistry*, 2nd ed., Cambridge University Press: Cambridge, 1984.

Lowry, T. H.; Richardson, K. S. *Mechanism and Theory in Organic Chemistry*, 3rd ed., Harper & Row: New York, 1987.

Moore, J. W.; Pearson, R. G. *Kinetics and Mechanism*, 3rd ed., John Wiley & Sons: New York, 1981.

Smith, M. B.; March, J. *Advanced Organic Chemistry*, 7th ed., John Wiley & Sons: New York, 2013.

These books include extensive bibliographies that permit the reader to delve more deeply into the subject.

Most libraries also subscribe to annual series of publications that specialize in articles dealing with reaction mechanisms. Among these are

Advances in Physical Organic Chemistry. Academic Press: London, 1963–present.

Annual Reports of the Chemical Society. Section B. Chemical Society: London, 1905–present. Specifically, the section Reaction Mechanisms.

Organic Reaction Mechanisms. John Wiley & Sons: Chichester, 1965–present.

Progress in Physical Organic Chemistry. Interscience: New York, 1963–present.

These publications provide the reader with citations from the original literature that can be very useful in an extensive literature search.

29.9 Organic Qualitative Analysis

Many laboratory manuals provide basic procedures for identifying organic compounds through a series of chemical tests and reactions. Occasionally, you might require a more complete description of analytical methods or a more complete set of tables of derivatives. Textbooks specializing in organic qualitative analysis should fill this need. Examples of sources for such information include

Cheronis, N. D.; Enriken, J. B. *Identification of Organic Compounds: A Student's Text Using Semimicro Techniques*. Interscience: New York, 1963.

Pasto, D. J.; Johnson, C. R. *Laboratory Text for Organic Chemistry: A Source Book of Chemical and Physical Techniques*. Prentice Hall: Englewood Cliffs, NJ, 1979.

Rappoport, Z. ed. *Handbook of Tables for Organic Compound Identification*, 3rd ed., CRC Press: Boca Raton, FL, 1967.

Shriner, R. L.; Hermann, C. K. F.; Merrill, T. C.; Curtin, D. Y.; Fuson, R. C. *The Systematic Identification of Organic Compounds*, 7th ed., John Wiley & Sons: New York, 1998.

Vogel, A. I. *Elementary Practical Organic Chemistry. Part 2. Qualitative Organic Analysis*, 2nd ed., John Wiley & Sons: New York, 1966.

Vogel, A. I. *Vogel's Textbook of Practical Organic Chemistry, Including Qualitative Organic Analysis*, 5th ed., Longman Group: London, 1989. Revised by members of the School of Chemistry, Thames Polytechnic.

29.10 Beilstein and Chemical Abstracts

One of the most useful sources of information about the physical properties, synthesis, and reactions of organic compounds is *Beilsteins Handbuch der Organischen Chemie*. This is a monumental work, initially edited by Friedrich Konrad Beilstein and updated through several revisions by the Beilstein Institute in Frankfurt am Main, Germany. The original edition (the *Hauptwerk*, abbreviated H) was published in 1918 and covers completely the literature to 1909. Five supplementary series (*Ergänzungswerken*) have been published since that time. The first supplement (*Erstes Ergänzungswerk*, abbreviated E I) covers the literature from 1910 to 1919; the second supplement (*Zweites Ergänzungswerk*, E II) covers 1920–1929; the third supplement (*Drittes Ergänzungswerk*, E III) covers 1930–1949; the fourth supplement (*Viertes Ergänzungswerk*, E IV) covers 1950–1959; and the fifth supplement (in English) covers 1960–1979. Volumes 17–27 of supplementary series III and IV, covering heterocyclic compounds, are combined in a joint issue, E III/IV. Supplementary series III, IV, and V are not complete, so the coverage of *Handbuch der Organischen Chemie* can be considered complete to 1929, with partial coverage to 1979.

Beilsteins Handbuch der Organischen Chemie, usually referred to simply as *Beilstein*, also contains two types of cumulative indices. The first of these is a name index (*Sachregister*), and the second is a formula index (*Formelregister*). These indices are particularly useful for a person wishing to locate a compound in *Beilstein*.

The principal difficulty in using *Beilstein* is that it is written in German through the fourth supplement. The fifth supplement is in English. Although some reading knowledge of German is useful, you can obtain information from the work by learning a few key phrases. For example, *Bildung* is “formation” or “structure.” *Darst* or *Darstellung* is “preparation,” *Kp* or *Siedepunkt* is “boiling point,” and *F* or *Schmelzpunkt* is “melting point.” Furthermore, the names of some compounds in German are not cognates of the English names. Some examples are *Apfelsäure* for “malic acid” (*säure* means “acid”), *Harnstoff* for “urea,” *Jod* for “iodine,” and *Zimtsäure* for “cinnamic acid.” If you have access to a German–English dictionary for chemists, many of these difficulties can be overcome. The best such dictionary is Patterson, A. M. *German–English Dictionary for Chemists*, 4th ed., John Wiley & Sons: New York, 1991.

Beilstein is organized according to a very sophisticated and complicated system. However, most students do not wish to become experts on *Beilstein* to this extent. A simpler, though slightly less reliable, method is to look for the compound in the formula index that accompanies the second supplement. By looking under the molecular formula, you will find the names of compounds that have that formula. After that name will be a series of numbers that indicate the pages and volume in which that compound is listed. Suppose, as an example, that you are searching for information on *p*-nitroaniline. This compound has the molecular formula $C_6H_6N_2O_2$. Searching for this formula in the formula index to the second supplement, you find

4-Nitro-anilin 12 711, I 349, II 383

This information tells you that *p*-nitroaniline is listed in the main edition, *Hauptwerk*, in Volume 12, page 711. Locate this particular volume, which is devoted to isocyclic monoamines and turn to page 711 to find the beginning of the section

on *p*-nitroaniline. At the left side of the top of this page is "Syst. No. 1671." This is the system number given to compounds in this part of Volume 12. The system number is useful, as it can help you find entries for this compound in subsequent supplements. The organization of *Beilstein* is such that all entries on *p*-nitroaniline in each of the supplements will be found in Volume 12. The entry in the formula index also indicates that material on this compound may be found in the first supplement on page 349 and in the second supplement on page 383. On page 349 of Volume 12 of the first supplement, there is a heading, "XII, 710–712," and on the left is "Syst. No. 1671." Material on *p*-nitroaniline is found in each supplement on a page that is headed with the volume and page of the *Hauptwerk* in which the same compound is found. On page 383 of Volume 12 of the second supplement, the heading in the center of the top of the page is "H12, 710–712." On the left, you find "Syst. No. 1671." Again, because *p*-nitroaniline appeared in Volume 12, page 711, of the main edition, you can locate it by searching through Volume 12 of any supplement until you find a page with the heading corresponding to Volume 12, page 711.

Because the third and fourth supplements are not complete, there is no comprehensive formula index for these supplements. However, you can still find material on *p*-nitroaniline by using the system number and the volume and page in the main work. In the third supplement, because the amount of information available has grown so much since the early days of *Beilstein's* work, Volume 12 has now expanded so that it is found in several bound parts. However, you select the part that includes system number 1671. In this part of Volume 12, you look through the pages until you find a page headed "Syst. No. 1671/H711." The information on *p*-nitroaniline is found on this page (page 1580). If Volume 12 of the fourth supplement were available, you would go on in the same way to locate more recent data on *p*-nitroaniline. This example is meant to illustrate how you can locate information on particular compounds without having to learn the *Beilstein* system of classification. You might do well to test your ability at finding compounds in *Beilstein* as we have described here.

Guidebooks to using *Beilstein*, which include a description of the *Beilstein* system, are recommended for anyone who wants to work extensively with *Beilstein*. Among such sources are:

Heller, S. R. *The Beilstein System: Strategies for Effective Searching*. Oxford University Press: New York, 1997.

How to Use Beilstein. Beilstein Institute, Frankfurt am Main. Springer-Verlag: Berlin, 1977.

Huntress, E. H. *A Brief Introduction to the Use of Beilsteins Handbuch der Organischen Chemie*, 2nd ed., John Wiley & Sons: New York, 1938.

Weissbach, O. *The Beilstein Guide: A Manual for the Use of Beilsteins Handbuch der Organischen Chemie*. Springer-Verlag: New York, 1976.

Beilstein reference numbers are listed in such handbooks as *CRC Handbook of Chemistry and Physics* and *Lange's Handbook of Chemistry*. Additionally, *Beilstein* numbers are included in the *Aldrich Handbook of Fine Chemicals*, issued by the Aldrich Chemical Company. If the compound you are seeking is listed in one of these handbooks, you will find that using *Beilstein* is simplified.

Another very useful publication for finding references for research on a particular topic is *Chemical Abstracts*, published by the Chemical Abstracts Service of the American Chemical Society. *Chemical Abstracts* contains abstracts of articles appearing in more than 10,000 journals from virtually every country conducting scientific research. These abstracts list the authors, the journal in which the article appeared,

the title of the article, and a short summary of the contents of the article. Abstracts of articles that appeared originally in a foreign language are provided in English, with a notation indicating the original language.

To use *Chemical Abstracts*, you must know how to use the various indices that accompany it. At the end of each volume, there appears a set of indices, including a formula index, a general subject index, a chemical substances index, an author index, and a patent index. The listings in each index refer the reader to the appropriate abstract according to the number assigned to it. There are also collective indices that combine all the indexed material appearing in a 5-year period (10-year period before 1956). In the collective indices, the listings include the volume number as well as the abstract number.

For material after 1929, *Chemical Abstracts* provides the most complete coverage of the literature. For material before 1929, use *Beilstein* before consulting *Chemical Abstracts*. *Chemical Abstracts* has the advantage that it is written entirely in English. Nevertheless, most students perform a literature search to find a relatively simple compound. Finding the desired entry for a simple compound is much easier in *Beilstein* than in *Chemical Abstracts*. For simple compounds, the indices in *Chemical Abstracts* are likely to contain very many entries. To locate the desired information, you must comb through this multitude of listings—potentially a very time-consuming task.

The opening pages of each index in *Chemical Abstracts* contain a brief set of instructions on using that index. If you want a more complete guide to *Chemical Abstracts*, consult a textbook designed to familiarize you with these abstracts and indices. Two such books are:

CAS Printed Access Tools: A Workbook. Chemical Abstracts Service, American Chemical Society: Washington, DC, 1977.

How to Search Printed CA. Chemical Abstracts Service, American Chemical Society: Washington, DC, 1989.

Chemical Abstracts Service maintains a computerized database that permits users to search through *Chemical Abstracts* rapidly and thoroughly. This service, which is called *CA Online*, is described in Section 29.11. *Beilstein* is also available for online searching by computer.

29.11 Computer Online Searching

You can search a number of chemistry databases online by using a computer and modem or a direct Internet connection. Many academic and industrial libraries can access these databases through their computers. One organization that maintains a large number of databases is the Scientific and Technical Information Network (STN International). The fee charged to the library for this service depends on the total time used in making the search, the type of information being asked for, the time of day when the search is being conducted, and the type of database being searched.

The Chemical Abstracts Service database (*CA Online*) is one of many databases available on STN. It is particularly useful to chemists. Unfortunately, this database extends back only to about 1967, although some earlier references are available. Searches for references earlier than 1967 must be made with printed abstracts (see Section 29.10). Searching online is much faster than searching in the printed abstracts. In addition, you can tailor the search in a number of ways by using keywords and the Chemical Abstracts Service Registry Number (CAS Number) as part of the search routine. The CAS Number is a specific number assigned to every compound listed in the *Chemical Abstracts* database. The CAS Number is used as a key in an online search to locate information about the compound. For the more

common organic compounds, you can easily obtain CAS Numbers from the catalogs of most of the companies that supply chemicals. Another advantage of performing an online search is that the *Chemical Abstracts* files are updated much more quickly than the printed versions of abstracts. This means that your search is more likely to reveal the most current information available.

Other useful databases available from STN include *Beilstein* and *CASREACTS*. As described in Section 29.10, *Beilstein* is very useful to organic chemists. Currently, there are more than 3.5 million compounds listed in the database. You can use the CAS Numbers to help in a search that has the potential of going back to 1830. *CASREACTS* is a chemical reactions database derived from over 100 journals covered by *Chemical Abstracts*, starting in 1985. With this database, you can specify a starting material and a product using the CAS Numbers. Further information on *CA Online*, *Beilstein*, *CASREACTS*, and other databases can be obtained from the following references:

Heller, S. R., ed. *The Beilstein Online Database: Implementation, Content and Retrieval*. American Chemical Society: Washington, DC, 1990.

Smith, M. B.; March, J. *Advanced Organic Chemistry*, 5th ed., John Wiley & Sons: New York, 2001.

Somerville, A. N. Information Sources for Organic Chemistry, 2: Searching by Functional Group. *Journal of Chemical Education*, 1991, 842.

Somerville, A. N. Subject Searching of Chemical Abstracts Online. *Journal of Chemical Education*, 1993, 200.

Wiggins, G. *Chemical Information Sources*. McGraw-Hill: New York, 1990. Integrates printed materials and computer sources of information.

SciFinder and SciFinder Scholar

The newest tools for online searching are SciFinder and SciFinder Scholar, the latter being the academic version of the software. This online service requires a yearly subscription and is available for use at many colleges and universities. SciFinder allows you to search several multidisciplinary CAS databases that contain information from as far back as 1907 to the present. The database may be searched in a variety of ways: by name, chemical substance, reaction, research topic, CAS number, or author. The program has drawing tools similar to ChemDraw, and what makes the program extremely useful is the ability to draw a structure on the screen and search for it. This avoids the need to name the structure first. In addition, substructure searching is allowed, which means that you may enter a partial structure and the program will find all references having compounds with the features you have indicated. Once you retrieve literature references, hyperlinks allow you view abstracts of the papers or retrieve physical property information. SciFinder is easy to use and requires minimal training. A recent book explains the program thoroughly:

Ridley, D. D. *Information Retrieval: SciFinder and SciFinder Scholar*. Wiley: New York, 2002.

For those who are at a university that subscribes to the service, there is an online tutorial at www.cas.org/SCIFINDER/SCHOLAR.

29.12 Scientific Journals

Ultimately, someone wanting information about a particular area of research will be required to read articles from the scientific journals. These journals are of two basic types: review journals and primary scientific journals. Journals that specialize in review articles summarize all of the work that bears on the particular topic. These articles may focus on the contributions of one particular researcher, but often

consider the contributions of many researchers to the subject. These articles also contain extensive bibliographies, which refer you to the original research articles. Among the important journals devoted, at least partly, to review articles are

Accounts of Chemical Research
Angewandte Chemie (International Edition, in English)
Chemical Reviews
Chemical Society Reviews (formerly known as *Quarterly Reviews*)
Nature
Science

The details of the research of interest appear in the primary scientific journals. Although there are thousands of journals published in the world, a few important journals specializing in articles dealing with organic chemistry include

Canadian Journal of Chemistry
European Journal of Organic Chemistry (formerly known as *Chemische Berichte*)
Journal of Organic Chemistry
Journal of the American Chemical Society
Journal of the Chemical Society, Chemical Communications
Journal of the Chemical Society, Perkin Transactions (Parts I and II)
Journal of Organometallic Chemistry
Organic Letters
Organometallics
Synlett
Synthesis
Tetrahedron
Tetrahedron Letters

29.13 Topics of Current Interest

The following journals and magazines are good sources for topics of educational and current interest. They specialize in news articles and focus on current events in chemistry or in science in general. Articles in these journals (magazines) can be useful in keeping you abreast of developments in science that are not part of your normal specialized scientific reading.

American Scientist
Chemical and Engineering News
Chemistry and Industry
Chemistry in Britain
Chemtech
Discover
Journal of Chemical Education
Nature
Omni
Science
Scientific American

Other sources for topics of current interest include the following:

Encyclopedia of Chemical Technology, 4th ed., 25 vols. plus index and supplements, 1992.

Also called *Kirk-Othmer Encyclopedia of Chemical Technology*.

McGraw-Hill Encyclopedia of Science and Technology, 20 volumes and supplements, 1997.

29.14 How to Conduct a Literature Search

The easiest method to follow in searching the literature is to begin with secondary sources and then go to the primary sources. In other words, you would try to locate material in a textbook, *Beilstein*, or *Chemical Abstracts*. From the results of that search, you would then consult one of the primary scientific journals.

A literature search that ultimately requires you to read one or more papers in the scientific journals is best conducted if you can identify a particular paper central to the study. Often, you can obtain this reference from a textbook or a review article on the subject. If this is not available, a search through *Beilstein* is required. A search through one of the handbooks that provides *Beilstein* reference numbers (see Section 29.10) may be helpful. Searching through *Chemical Abstracts* would be considered the next logical step. From these sources, you should be able to identify citations from the original literature on the subject.

Additional citations may be found in the references cited in the journal article. In this way, the background leading to the research can be examined. It is also possible to conduct a search forward in time from the date of the journal article through the *Science Citation Index*. This publication provides the service of listing articles and the papers in which these articles were cited. Although the *Science Citation Index* consists of several types of indices, the *Citation Index* is most useful for the purposes described here. A person who knows of a particular key reference on a subject can examine the *Science Citation Index* to obtain a list of papers that have used that seminal reference in support of the work described. The *Citation Index* lists papers by their senior author, journal, volume, page, and date, followed by citations of papers that have referred to that article, author, journal, volume, page, and date of each. The *Citation Index* is published in annual volumes, with quarterly supplements issued during the current year. Each volume contains a complete list of the citations of the key articles made during that year. A disadvantage is that *Science Citation Index* has been available only since 1961. An additional disadvantage is that you may miss journal articles on the subject of interest if *Citation Index* failed to cite that particular key reference in its bibliographies—a reasonably likely possibility.

You can, of course, conduct a literature search by a “brute force” method, by beginning the search with *Beilstein* or even with the indices in *Chemical Abstracts*. However, the task can be made much easier by performing a computer search (see Section 29.11) or by starting with a book or an article of general and broad coverage, which can provide a few citations for starting points in the search.

The following guides to using the chemical literature are provided for the reader who is interested in going farther into this subject.

Bottle, R. T.; Rowland, J. F. B., eds. *Information Sources in Chemistry*, 4th ed., Bowker-Saur: New York, 1992.

Maizell, R. E. *How to Find Chemical Information: A Guide for Practicing Chemists, Educators, and Students*, 3rd ed., John Wiley & Sons: New York, 1998.

Mellon, M. G. *Chemical Publications*, 5th ed., McGraw-Hill: New York, 1982.

Wiggins, G. *Chemical Information Sources*. McGraw-Hill: New York, 1991. Integrates printed materials and computer sources of information.

PROBLEMS

1. Find the following compounds in the formula index for the *Second Supplement of Beilstein* (see Section 29.10). (1) List the page numbers from the main work and the supplements (first and second). (2) Using these page numbers, look up the

- system number (Syst. No.) and the main work number (*Hauptwerk* number, H) for each compound in the main work and the first and second supplements. In some cases, a compound may not be found in all three places. (3) Now use the system number and main work number to find each of these compounds in the third and fourth supplements. List the page numbers where these compounds are found.
- 2,5-hexanedione (acetonylacetone)
 - 3-nitroacetophenone
 - 4-*tert*-butylcyclohexanone
 - 4-phenylbutanoic acid (4-phenylbutyric acid, γ -phenylbuttersäure)
- Using the *Science Citation Index* (see Section 29.14), list five research papers by complete title and journal citation for each of the following chemists who have been awarded the Nobel Prize. Use the *Five-Year Cumulative Source Index* for the years 1980–1984 as your source.
 - H. C. Brown
 - R. B. Woodward
 - D. J. Cram
 - G. Olah
 - The reference book by Smith and March is listed in Section 29.2. Using Appendix 2 in this book, give two methods for preparing the following functional groups. You will need to provide equations.
 - carboxylic acids
 - aldehydes
 - esters (carboxylic esters)
 - Organic Syntheses* is described in Section 29.6. There are currently nine collective volumes in the series, each with its own index. Find the compounds listed below and provide the equations for preparing each compound.
 - 2-methylcyclopentane-1,3-dione
 - cis*- Δ^4 -tetrahydrophthalic anhydride (listed as tetrahydrophthalic anhydride)
 - Provide four methods that may be used to oxidize an alcohol to an aldehyde. Give complete literature references for each method, as well as equations. Use the *Compendium of Organic Synthetic Methods* or *Survey of Organic Syntheses* by Buehler and Pearson (see Section 29.6).

Appendices

Tables of Unknowns and Derivatives

More extensive tables of unknowns may be found in Z. Rappoport, ed. *Handbook of Tables for Organic Compound Identification*, 3rd ed. CRC Press: Boca Raton FL, 1967.

ALDEHYDES

Compound	BP	MP	Semi-carbazone*	2,4-Dinitro-phenyl-hydrazone*
Ethanal (acetaldehyde)	21	—	162	168
Propanal (propionaldehyde)	48	—	89	148
Propenal (acrolein)	52	—	171	165
2-Methylpropanal (isobutyraldehyde)	64	—	125	187
Butanal (butyraldehyde)	75	—	95	123
3-Methylbutanal (isovaleraldehyde)	92	—	107	123
Pentanal (valeraldehyde)	102	—	—	106
2-Butenal (crotonaldehyde)	104	—	199	190
2-Ethylbutanal (diethylacetaldehyde)	117	—	99	95
Hexanal (caproaldehyde)	130	—	106	104
Heptanal (heptaldehyde)	153	—	109	108
2-Furaldehyde (furfural)	162	—	202	212
2-Ethylhexanal	163	—	254	114
Octanal (caprylaldehyde)	171	—	101	106
Benzaldehyde	179	—	222	237
Nonanal (nonyl aldehyde)	185	—	100	100
Phenylethanal (phenylacetaldehyde)	195	33	153	121
2-Hydroxybenzaldehyde (salicylaldehyde)	197	—	231	248
4-Methylbenzaldehyde (<i>p</i> -tolualdehyde)	204	—	234	234
3,7-Dimethyl-6-octenal (citronellal)	207	—	82	77
Decanal (decyl aldehyde)	207	—	102	104
2-Chlorobenzaldehyde	213	11	229	213
3-Chlorobenzaldehyde	214	18	228	248
3-Methoxybenzaldehyde (<i>m</i> -anisaldehyde)	230	—	233 d.	—
3-Bromobenzaldehyde	235	—	205	—
4-Methoxybenzaldehyde (<i>p</i> -anisaldehyde)	248	2.5	210	253
<i>trans</i> -Cinnamaldehyde	250 d.	—	215	255

ALDEHYDES (Cont.)

Compound	BP	MP	Semi-carbazone*	2,4-Dinitro-phenyl-hydrazone*
3,4-Methylenedioxybenzaldehyde (piperonal)	263	37	230	266 d.
2-Methoxybenzaldehyde (<i>o</i> -anisaldehyde)	245	38	215 d.	254
3,4-Dimethoxybenzaldehyde	—	44	177	261
2-Nitrobenzaldehyde	—	44	256	265
4-Chlorobenzaldehyde	—	48	230	254
4-Bromobenzaldehyde	—	57	228	257
3-Nitrobenzaldehyde	—	58	246	293
2,4-Dimethoxybenzaldehyde	—	71	—	—
2,4-Dichlorobenzaldehyde	—	72	—	—
4-Dimethylaminobenzaldehyde	—	74	222	325
4-Hydroxy-3-methoxybenzaldehyde (vanillin)	—	82	230	271
3-Hydroxybenzaldehyde	—	104	198	259
5-Bromo-2-hydroxybenzaldehyde (5-bromosalicylaldehyde)	—	106	297 d.	—
4-Nitrobenzaldehyde	—	106	221	320 d.
4-Hydroxybenzaldehyde	—	116	224	280 d.
(±)-Glyceraldehyde	—	142	160 d.	167

Note: "d" indicates "decomposition."

*See Appendix 2, "Procedures for Preparing Derivatives."

KETONES

Compound	BP	MP	Semi-carbazone*	2,4-Dinitro-phenyl-hydrazone*
2-Propanone (acetone)	56	—	187	126
2-Butanone (methyl ethyl ketone)	80	—	146	117
3-Buten-2-one (methyl vinyl ketone)	81	—	140	—
3-Methyl-2-butanone (isopropyl methyl ketone)	94	—	112	120
2-Pentanone (methyl propyl ketone)	102	—	112	143
3-Pentanone (diethyl ketone)	102	—	138	156
3,3-Dimethyl-2-butanone (pinacolone)	106	—	157	125
4-Methyl-2-pentanone (isobutyl methyl ketone)	117	—	132	95
2,4-Dimethyl-3-pentanone (diisopropyl ketone)	124	—	160	86
3-Hexanone	125	—	113	130
2-Hexanone (methyl butyl ketone)	128	—	121	106
4-Methyl-3-penten-2-one (mesityl oxide)	130	—	164	200
Cyclopentanone	131	—	210	146

KETONES (Cont.)

Compound	BP	MP	Semi-carbazone*	2,4-Dinitro-phenyl-hydrazone*
5-Hexen-2-one	131	—	102	108
2,3-Pentanedione	134	—	122 (mono) 209 (di)	209
5-Methyl-3-hexanone	136	—	—	—
2,4-Pentanedione (acetylacetone)	139	—	122 (mono) 209 (di)	209
4-Heptanone (dipropyl ketone)	144	—	132	75
5-Methyl-2-hexanone	145	—	—	—
1-Hydroxy-2-propanone (hydroxyacetone, acetol)	146	—	196	129
3-Heptanone	148	—	101	—
2-Heptanone (methyl amyl ketone)	151	—	123	89
Cyclohexanone	156	—	166	162
2-Methylcyclohexanone	165	—	191	136
3-Octanone	167	—	—	—
2,6-Dimethyl-4-heptanone (diisobutyl ketone)	168	—	122	66
2-Octanone	173	—	122	92
Cycloheptanone	181	—	163	148
Ethyl acetoacetate	181	—	129 d.	93
5-Nonanone	186	—	90	—
3-Nonanone	187	—	112	—
2,5-Hexanedione (acetylacetone)	191	-9	185 (mono) 244 (di)	257 (di)
2-Nonanone	195	-8	118	—
Acetophenone (methyl phenyl ketone)	202	20	198	238
2-Hydroxyacetophenone	215	28	210	212
1-Phenyl-2-propanone (phenylacetone)	216	27	198	156
Propiophenone (1-phenyl-1-propanone)	218	21	173	191
Isobutyrophenone (2-methyl-1-phenyl-1-propanone)	222	—	181	163
1-Phenyl-2-butanone	226	—	135	—
4-Methylacetophenone	226	28	205	258
3-Chloroacetophenone	228	—	232	—
2-Chloroacetophenone	229	—	160	—
Butyrophenone (1-phenyl-1-butanone)	230	12	187	190
2-Undecanone	231	12	122	63
4-Chloroacetophenone	232	12	204	231
4-Phenyl-2-butanone (benzylacetone)	235	—	142	127
2-Methoxyacetophenone	239	—	183	—
3-Methoxyacetophenone	240	—	196	—
Valerophenone (1-phenyl-1-pentanone)	248	—	160	166

KETONES (Cont.)

Compound	BP	MP	Semi-carbazone*	2,4-Dinitrophenylhydrazone*
4-Chloropropiophenone	—	36	176	—
4-Phenyl-3-buten-2-one (benzalacetone)	—	37	187	227
4-Methoxyacetophenone	—	38	198	220
3-Bromopropiophenone	—	40	183	—
1-Indanone	—	41	233	258
Benzophenone	—	48	164	238
4-Bromoacetophenone	—	51	208	230
3,4-Dimethoxyacetophenone	—	51	218	207
2-Acetonaphthone (methyl 2-naphthyl ketone)	—	53	234	262 d.
Desoxybenzoin (benzyl phenyl ketone)	—	60	148	204
1,1-Diphenylacetone	—	61	170	—
4-Chlorobenzophenone	—	76	—	185
3-Nitroacetophenone	—	80	257	228
4-Nitroacetophenone	—	80	—	—
4-Bromobenzophenone	—	82	350	230
Fluorenone	—	83	—	283
4-Hydroxyacetophenone	—	109	199	210
Benzoin	—	136	206	245
4-Hydroxypropiophenone	—	148	—	229
(±)-Camphor	—	179	237	164

Note: "d" indicates "decomposition."

*See Appendix 2, "Procedures for Preparing Derivatives."

CARBOXYLIC ACIDS

Compound	BP	MP	<i>p</i> -Toluidide*	Anilide*	Amide*
Methanoic acid (formic acid)	101	8	53	47	43
Ethanoic acid (acetic acid)	118	17	148	114	82
Propenoic acid (acrylic acid)	139	13	141	104	85
Propanoic acid (propionic acid)	141	—	124	103	81
2-Methylpropanoic acid (isobutyric acid)	154	—	104	105	128
Butanoic acid (butyric acid)	162	—	72	95	115
3-Butenoic acid (vinylacetic acid)	163	—	—	58	73
2-Methylpropenoic acid (methacrylic acid)	163	16	—	87	102
Pyruvic acid	165 d.	14	109	104	124

CARBOXYLIC ACIDS (Cont.)

Compound	BP	MP	<i>p</i> -Toluidide*	Anilide*	Amide*
3-Methylbutanoic acid (isovaleric acid)	176	—	106	109	135
3,3-Dimethylbutanoic acid	185	—	134	132	132
Pentanoic acid (valeric acid)	186	—	74	63	106
2-Chloropropanoic acid	186	—	124	92	80
Dichloroacetic acid	194	6	153	118	98
2-Methylpentanoic acid	195	—	80	95	79
Hexanoic acid (caproic acid)	205	—	75	95	101
2-Bromopropanoic acid	205 d.	24	125	99	123
Heptanoic acid	223	—	81	70	96
2-Ethylhexanoic acid	228	—	—	—	102
Cyclohexanecarboxylic acid	233	31	—	146	186
Octanoic acid (caprylic acid)	237	16	70	57	107
Nonanoic acid	254	12	84	57	99
Decanoic acid (capric acid)	—	32	78	70	108
4-Oxopentanoic acid (levulinic acid)	—	33	108	102	108 d.
Trimethylacetic acid (pivalic acid)	—	35	120	130	155
3-Chloropropanoic acid	—	40	—	—	101
Dodecanoic acid (lauric acid)	—	43	87	78	100
3-Phenylpropanoic acid (hydrocinnamic acid)	—	48	135	98	105
Bromoacetic acid	—	50	—	131	91
4-Phenylbutanoic acid	—	52	—	—	84
Tetradecanoic acid (myristic acid)	—	54	93	84	103
Trichloroacetic acid	—	57	113	97	141
3-Bromopropanoic acid	—	61	—	—	111
Hexadecanoic acid (palmitic acid)	—	62	98	90	106
Chloroacetic acid	—	63	162	137	121
Cyanoacetic acid	—	66	—	198	120
Octadecanoic acid (stearic acid)	—	69	102	95	109
<i>trans</i> -2-Butenoic acid (crotonic acid)	—	72	132	118	158
Phenylacetic acid	—	77	136	118	156
α -Methyl- <i>trans</i> -cinnamic acid	—	81	—	—	128
4-Methoxyphenylacetic acid	—	87	—	—	189
3,4-Dimethoxyphenyl acetic acid	—	97	—	—	147
Pentanedioic acid (glutaric acid)	—	98	218 (di)	224 (di)	176 (di)
Phenoxyacetic acid	—	99	—	99	102
2-Methoxybenzoic acid (<i>o</i> -anisic acid)	—	100	—	131	129
2-Methylbenzoic acid (<i>o</i> -toluic acid)	—	104	144	125	142
Nonanedioic acid (azelaic acid)	—	106	201 (di)	107 (mono) 186 (di)	93 (mono) 175 (di)

CARBOXYLIC ACIDS (Cont.)

Compound	BP	MP	<i>p</i> -Toluidide*	Anilide*	Amide*
3-Methoxybenzoic acid (<i>m</i> -anisic acid)	—	107	—	—	136
3-Methylbenzoic acid (<i>m</i> -toluic acid)	—	111	118	126	94
4-Bromophenylacetic acid	—	117	—	—	194
(±)-Phenylhydroxyacetic acid (mandelic acid)	—	118	172	151	133
Benzoic acid	—	122	158	163	130
2,4-Dimethylbenzoic acid	—	126	—	141	180
2-Benzoylbenzoic acid	—	127	—	195	165
Maleic acid	—	130	142 (di)	198 (mono) 187 (di)	172 (mono) 260 (di)
Decanedioic acid (sebacic acid)	—	133	201 (di)	122 (mono) 200 (di)	170 (mono) 210 (di)
3-Chlorocinnamic acid	—	133	142	135	76
2-Furoic acid	—	133	170	124	143
<i>trans</i> -Cinnamic acid	—	133	168	153	147
2-Acetylsalicylic acid (aspirin)	—	138	—	136	138
5-Chloro-2-nitrobenzoic acid	—	139	—	164	154
2-Chlorobenzoic acid	—	140	131	118	139
3-Nitrobenzoic acid	—	140	162	155	143
4-Chloro-2-nitrobenzoic acid	—	142	—	—	172
2-Nitrobenzoic acid	—	146	—	155	176
2-Aminobenzoic acid (anthranilic acid)	—	146	151	131	109
Diphenylacetic acid	—	148	172	180	167
2-Bromobenzoic acid	—	150	—	141	155
Benzilic acid	—	150	190	175	154
Hexanedioic acid (adipic acid)	—	152	239	151 (mono) 241 (di)	125 (mono) 220 (di)
Citric acid	—	153	189 (tri)	198 (tri)	210 (tri)
4-Nitrophenylacetic acid	—	153	—	198	198
2,5-Dichlorobenzoic acid	—	153	—	—	155
3-Chlorobenzoic acid	—	156	—	123	134
2,4-Dichlorobenzoic acid	—	158	—	—	194
4-Chlorophenoxyacetic acid	—	158	—	125	133
2-Hydroxybenzoic acid (salicylic acid)	—	158	156	136	142
5-Bromo-2-hydroxybenzoic acid (5-bromosalicylic acid)	—	165	—	222	232
3,4-Dimethylbenzoic acid	—	165	—	104	130
2-Chloro-5-nitrobenzoic acid	—	166	—	—	178
Methylenesuccinic acid (itaconic acid)	—	166 d.	—	152 (mono)	191 (di)
(+)-Tartaric acid	—	169	—	180 (mono) 264 (di)	171 (mono) 196 (di)

CARBOXYLIC ACIDS (Cont.)

Compound	BP	MP	<i>p</i> -Toluidide*	Anilide*	Amide*
5-Chlorosalicylic acid	—	172	—	—	227
4-Methylbenzoic acid (<i>p</i> -toluic acid)	—	180	160	145	160
4-Chloro-3-nitrobenzoic acid	—	182	—	131	156
4-Methoxybenzoic acid (<i>p</i> -anisic acid)	—	184	186	169	167
Butanedioic acid (succinic acid)	—	188	180 (mono) 255 (di)	143 (mono) 230 (di)	157 (mono) 260 (di)
4-Ethoxybenzoic acid	—	198	—	170	202
Fumaric acid	—	200 s.	—	233 (mono) 314 (di)	270 (mono) 266 (di)
3-Hydroxybenzoic acid	—	201 s.	163	157	170
3,5-Dinitrobenzoic acid	—	202	—	234	183
3,4-Dichlorobenzoic acid	—	209	—	—	133
Phthalic acid	—	210 d.	150 (mono) 201 (di)	169 (mono) 253 (di)	144 (mono) 220 (di)
4-Hydroxybenzoic acid	—	214	204	197	162
3-Nitrophthalic acid	—	215	226 (di)	234 (di)	201 (di)
Pyridine-3-carboxylic acid (nicotinic acid)	—	236	150	132	128
4-Nitrobenzoic acid	—	240	204	211	201
4-Chlorobenzoic acid	—	242	—	194	179
4-Bromobenzoic acid	—	251	—	197	190

Note: "d" indicates "decomposition"; "s" indicates "sublimation."

*See Appendix 2, "Procedures for Preparing Derivatives."

PHENOLS†

Compound	BP	MP	α -Naphthyl- urethane*	Bromo Derivative*			
				Mono	Di	Tri	Tetra
2-Chlorophenol	176	7	120	48	76	—	—
3-Methylphenol (<i>m</i> -cresol)	203	12	128	—	—	84	—
2-Ethylphenol	207	—	—	—	—	—	—
2,4-Dimethylphenol	212	23	135	—	—	—	—
2-Methylphenol (<i>o</i> -cresol)	191	32	142	—	56	—	—
2-Methoxyphenol (guaiacol)	204	32	118	—	—	116	—
4-Methylphenol (<i>p</i> -cresol)	202	35	146	—	49	—	198
3-Chlorophenol	214	35	158	—	—	—	—
4-Methyl-2-nitrophenol	—	35	—	—	—	—	—
2,4-Dibromophenol	238	40	—	95	—	—	—
Phenol	181	42	133	—	—	95	—
4-Chlorophenol	217	43	166	33	90	—	—

PHENOLS[†] (Cont.)

Compound	BP	MP	α -Naphthyl-urethane*	Bromo Derivative*			
				Mono	Di	Tri	Tetra
4-Ethylphenol	219	45	128	—	—	—	—
2-Nitrophenol	216	45	113	—	117	—	—
2-Isopropyl-5-methylphenol (thymol)	234	51	160	55	—	—	—
4-Methoxyphenol	243	56	—	—	—	—	—
3,4-Dimethylphenol	225	64	141	—	—	171	—
4-Bromophenol	238	64	169	—	—	—	—
4-Chloro-3-methylphenol	235	66	153	—	—	—	—
3,5-Dimethylphenol	220	68	—	—	—	166	—
2,6-Di- <i>tert</i> -butyl-4-methylphenol	—	70	—	—	—	—	—
2,4,6-Trimethylphenol	232	72	—	—	—	—	—
2,5-Dimethylphenol	212	75	173	—	—	178	—
1-Naphthol (α -naphthol)	278	94	152	—	105	—	—
2-Methyl-4-nitrophenol	186	96	—	—	—	—	—
2-Hydroxyphenol (catechol)	245	104	175	—	—	—	192
2-Chloro-4-nitrophenol	—	106	—	—	—	—	—
3-Hydroxyphenol (resorcinol)	—	109	—	—	—	112	—
4-Nitrophenol	—	112	150	—	142	—	—
2-Naphthol (β -naphthol)	—	123	157	84	—	—	—
3-Methyl-4-nitrophenol	—	129	—	—	—	—	—
1,2,3-Trihydroxybenzene (pyrogallol)	—	133	—	—	158	—	—
4-Phenylphenol	—	164	—	—	—	—	—

*See Appendix 2, "Procedures for Preparing Derivatives."

[†]Also check:

Salicylic acid (2-hydroxybenzoic acid)

Esters of salicylic acid (salicylates)

Salicylaldehyde (2-hydroxybenzaldehyde)

4-Hydroxybenzaldehyde

4-Hydroxypropiophenone

3-Hydroxybenzoic acid

4-Hydroxybenzoic acid

4-Hydroxybenzophenone

PRIMARY AMINES[†]

Compound	BP	MP	Benzamide*	Acetamide*
<i>t</i> -Butylamine	46	—	134	101
Propylamine	48	—	84	—
Allylamine	56	—	—	—
<i>sec</i> -Butylamine	63	—	76	—

PRIMARY AMINES[†] (Cont.)

Compound	BP	MP	Benzamide*	Acetamide*
Isobutylamine	69	—	57	—
Butylamine	78	—	42	—
Isopentylamine (ioamylamine)	96	—	—	—
Pentylamine (amylamine)	104	—	—	—
Enthylenediamine	118	—	244 (di)	172 (di)
Hexylamine	132	—	40	—
Cyclohexylamine	135	—	149	101
1,3-Diaminopropane	140	—	148 (di)	126 (di)
Furfurylamine	145	—	—	—
Heptylamine	156	—	—	—
Octylamine	180	—	—	—
Benzylamine	184	—	105	65
Aniline	184	—	163	114
2-Methylaniline (<i>o</i> -toluidine)	200	—	144	110
3-Methylaniline (<i>m</i> -toluidine)	203	—	125	65
2-Chloroaniline	208	—	99	87
2,6-Dimethylaniline	216	11	168	177
2,5-Dimethylaniline	216	14	140	139
3,5-Dimethylaniline	220	—	144	—
4-Isopropylaniline	225	—	162	102
2-Methoxyaniline (<i>o</i> -anisidine)	225	6	60	85
3-Chloroaniline	230	—	120	74
2-Ethoxyaniline (<i>o</i> -phenetidine)	231	—	104	79
4-Chloro-2-methylaniline	241	29	142	140
4-Ethoxyaniline (<i>p</i> -phenetidine)	250	2	173	137
3-Bromoaniline	251	18	120	87
2-Bromoaniline	250	31	116	99
2,6-Dichloroaniline	—	39	—	—
4-Methylaniline (<i>p</i> -toluidine)	200	43	158	147
2-Ethylaniline	210	47	147	111
2,5-Dichloroaniline	251	50	120	132
4-Methoxyaniline (<i>p</i> -anisidine)	—	58	154	130
2,4-Dichloroaniline	245	62	117	145
4-Bromoaniline	245	64	204	168
4-Chloroaniline	—	72	192	179
2-Nitroaniline	—	72	110	92
2,4,6-Trichloroaniline	262	75	174	204
Ethyl <i>p</i> -aminobenzoate	—	89	148	110
<i>o</i> -Phenylenediamine	258	102	301 (di)	185 (di)
2-Methyl-5-nitroaniline	—	106	186	151
4-Aminoacetophenone	—	106	205	167

PRIMARY AMINES[†] (Cont.)

Compound	BP	MP	Benzamide*	Acetamide*
2-Chloro-4-nitroaniline	—	108	161	139
3-Nitroaniline	—	114	157	155
4-Methyl-2-nitroaniline	—	116	148	99
4-Chloro-2-nitroaniline	—	118	133	104
2,4,6-Tribromoaniline	—	120	200	232
2-Methyl-4-nitroaniline	—	130	—	202
2-Methoxy-4-nitroaniline	—	138	149	153
<i>p</i> -Phenylenediamine	—	140	128 (mono) 300 (di)	162 (mono) 304 (di)
4-Nitroaniline	—	148	199	215
4-Aminoacetanilide	—	162	—	304
2,4-Dinitroaniline	—	180	202	120

*See Appendix 2, "Procedures for Preparing Derivatives."

†Also check 4-aminobenzoic acid and its esters.

SECONDARY AMINES

Compound	BP	MP	Benzamide*	Acetamide*
Diethylamine	56	—	42	—
Diisopropylamine	84	—	—	—
Pyrrolidine	88	—	Oil	—
Piperidine	106	—	48	—
Dipropylamine	110	—	—	—
Morpholine	129	—	75	—
Diisobutylamine	139	—	—	86
<i>N</i> -Methylcyclohexylamine	148	—	85	—
Dibutylamine	159	—	—	—
Benzylmethylamine	184	—	—	—
<i>N</i> -Methylaniline	196	—	63	102
<i>N</i> -Ethylaniline	205	—	60	54
<i>N</i> -Ethyl- <i>m</i> -toluidine	221	—	72	—
Dicyclohexylamine	256	—	153	103
<i>N</i> -Benzylaniline	298	37	107	58
Indole	254	52	68	157
Diphenylamine	302	52	180	101
<i>N</i> -phenyl-1-naphthylamine	335	62	152	115

*See Appendix 2, "Procedures for Preparing Derivatives."

TERTIARY AMINES†

Compound	BP	MP	Methiodide*
Triethylamine	89	—	280
Pyridine	115	—	117
2-Methylpyridine (α -picoline)	129	—	230
2,6-Dimethylpyridine (2,6-lutidine)	143	—	233
4-Methylpyridine (4-picoline)	143	—	—
3-Methylpyridine (β -picoline)	144	—	92
Tripopylamine	157	—	207
<i>N,N</i> -Dimethylbenzylamine	183	—	179
<i>N,N</i> -Dimethylaniline	193	—	228 d.
Tributylamine	216	—	186
<i>N,N</i> -Diethylaniline	217	—	102
Quinoline	237	—	72/133

Note: "d" indicates "decomposition."

*See Appendix 2, "Procedures for Preparing Derivatives."

† Also check nicotinic acid and its esters.

ALCOHOLS

Compound	BP	MP	3,5-Di-nitrobenzoate*	Phenyl-urethane*
Methanol	65	—	108	47
Ethanol	78	—	93	52
2-Propanol (isopropyl alcohol)	82	—	123	88
2-Methyl-2-propanol (<i>t</i> -butyl alcohol)	83	26	142	136
3-Buten-2-ol	96	—	54	—
2-Propen-1-ol (allyl alcohol)	97	—	49	70
1-Propanol	97	—	74	57
2-Butanol (<i>sec</i> -butyl alcohol)	99	—	76	65
2-Methyl-2-butanol (<i>t</i> -pentyl alcohol)	102	-8.5	116	42
2-Methyl-3-butyn-2-ol	104	—	112	—
2-Methyl-1-propanol (isobutyl alcohol)	108	—	87	86
3-Buten-1-ol	113	—	59	25
3-Methyl-2-butanol	114	—	76	68
2-Propyn-1-ol (propargyl alcohol)	114	—	—	—
3-Pentanol	115	—	101	48
1-Butanol	118	—	64	61
2-Pentanol	119	—	62	—
3,3-Dimethyl-2-butanol	120	—	107	77
2,3-Dimethyl-2-butanol	121	—	111	65
2-Methyl-2-pentanol	123	—	72	—

ALCOHOLS (Cont.)

Compound	BP	MP	3,5-Di-nitrobenzoate*	Phenyl-urethane*
3-Methyl-3-pentanol	123	—	96	43
2-Methoxyethanol	124	—	—	(113) [†]
2-Methyl-3-pentanol	128	—	85	50
2-Chloroethanol	129	—	95	51
3-Methyl-1-butanol (isoamyl alcohol)	132	—	61	56
4-Methyl-2-pentanol	132	—	65	143
2-Ethoxyethanol	135	—	75	(67) [†]
3-Hexanol	136	—	97	—
1-Pentanol	138	—	46	46
2-Hexanol	139	—	39	(61) [†]
2,4-Dimethyl-3-pentanol	140	—	—	95
Cyclopentanol	140	—	115	132
2-Ethyl-1-butanol	146	—	51	—
2,2,2-Trichloroethanol	151	—	142	87
1-Hexanol	157	—	58	42
2-Heptanol	159	—	49	(54) [†]
Cyclohexanol	160	—	113	82
3-Chloro-1-propanol	161	—	77	38
(2-Furyl)-methanol (furfuryl alcohol)	170	—	80	45
1-Heptanol	176	—	47	60
2-Octanol	179	—	32	114
2-Ethyl-1-hexanol	185	—	—	(61) [†]
1-Octanol	195	—	61	74
3,7-Dimethyl-1,6-octadien-3-ol (linalool)	196	—	—	66
2-Nonanol	198	—	43	(56) [†]
Benzyl alcohol	204	—	113	77
1-Phenylethanol	204	20	92	95
1-Nonanol	214	—	52	62
1,3-Propanediol	215	—	178 (di)	137 (di)
2-Phenylethanol	219	—	108	78
1-Decanol	231	7	57	59
3-Phenylpropanol	236	—	45	92
1-Dodecanol (lauryl alcohol)	—	24	60	74
3-Phenyl-2-propen-1-ol (cinnamyl alcohol)	250	34	121	90
α -Terpineol	221	36	78	112
1-Tetradecanol (myristyl alcohol)	—	39	67	74
(-)-Menthol	212	41	158	111
1-Hexadecanol (cetyl alcohol)	—	49	66	73
2,2-Dimethyl-1-propanol (neopentyl alcohol)	113	56	—	144
4-Methylbenzyl alcohol	217	59	117	79

ALCOHOLS (Cont.)

Compound	BP	MP	3,5-Di-nitrobenzoate*	Phenyl-urethane*
1-Octadecanol (stearyl alcohol)	—	59	77	79
Diphenylmethanol (benzhydrol)	—	68	141	139
4-Nitrobenzyl alcohol	—	93	157	—
Benzoin	—	136	—	165
Cholesterol	—	147	—	168
Triphenylmethanol	—	161	—	—
(+)-Borneol	—	208	154	138

*See Appendix 2, "Procedures for Preparing Derivatives."

† α -Naphthylurethane.**ESTERS**

Compound	BP	MP	Compound	BP	MP
Methyl formate	32	—	(ethyl 2-methylpropanoate)	110	—
Ethyl formate	54	—	Isopropyl propionate		
Methyl acetate	57	—	(isopropyl propanoate)	110	—
Isopropyl formate	71	—	2-Butyl acetate (<i>sec</i> -butyl acetate)	111	—
Vinyl acetate	72	—	Methyl isovalerate	117	—
Ethyl acetate	77	—	(methyl 3-methylbutanoate)		
Methyl propionate			Isobutyl acetate		
(methyl propanoate)	80	—	(2-methylpropyl acetate)	117	—
Methyl acrylate	80	—	Ethyl pivalate		
Propyl formate	81	—	(ethyl 2,2-dimethylpropanoate)	118	—
Isopropyl acetate	89	—	Methyl crotonate (methyl 2-butenate)	119	—
Ethyl chloroformate	93	—	Ethyl butyrate (ethyl butanoate)	121	—
Methyl isobutyrate	93	—	Propyl propionate (propyl propanoate)	123	—
(methyl 2-methylpropanoate)			Butyl acetate	126	—
2-Propenyl acetate (isopropenyl acetate)	94	—	Methyl valerate (methyl pentanoate)	128	—
<i>tert</i> -Butyl acetate			Methyl methoxyacetate	130	—
(1,1-dimethylethyl acetate)	98	—	Methyl chloroacetate	130	—
Ethyl propionate (ethyl propanoate)	99	—	Ethyl isovalerate		
Methyl methacrylate			(ethyl 3-methylbutanoate)	134	—
(methyl 2-methylpropenoate)	100	—	Ethyl crotonate (ethyl 2-butenate)	138	—
Methyl pivalate			Isopentyl acetate		
(methyl trimethyl acetate)	101	—	(3-methylbutyl acetate)	142	—
Ethyl acrylate (ethyl propenoate)	101	—	2-Methoxyethyl acetate	145	—
Propyl acetate	102	—	Ethyl chloroacetate	145	—
Methyl butyrate (methyl butanoate)	102	—	Ethyl valerate (ethyl pentanoate)	146	—
Ethyl isobutyrate			Ethyl α -chloropropanoate	146	—

ESTERS (Cont.)

Compound	BP	MP	Compound	BP	MP
Pentyl acetate	147	—	Methyl salicylate	224	—
Methyl hexanoate	151	—	Diethyl maleate	224	—
Ethyl lactate	154	—	Ethyl phenylacetate	228	—
Butyl butyrate	167	—	Propyl benzoate	231	—
Ethyl hexanoate	168	—	Ethyl salicylate	234	—
Hexyl acetate	169	—	Dimethyl suberate	268	—
Methyl acetoacetate	170	—	Ethyl cinnamate	271	—
Methyl heptanoate (methyl enanthlate)	172	—	Dimethyl phthalate	284	—
Furfuryl acetate	176	—	Diethyl phthalate	298	—
Methyl 2-furoate	181	—	Methyl cinnamate	—	36
Dimethyl malonate	181	—	Ethyl 2-furoate	—	36
Ethyl acetoacetate	181	—	Methyl stearate	—	39
Diethyl oxalate	185	—	Dimethyl itaconate	—	39
Ethyl heptanoate	187	—	Phenyl salicylate	—	42
Heptyl acetate	192	—	Diethyl terephthalate	—	44
Dimethyl succinate	196	—	Methyl 4-chlorobenzoate	—	44
Phenyl acetate	197	—	Ethyl 3-nitrobenzoate	—	47
Diethyl malonate	199	—	Methyl mandelate	—	53
Methyl benzoate	199	—	Ethyl 4-nitrobenzoate	—	56
Dimethyl maleate	204	—	Dimethyl isophthalate	—	68
Ethyl levulinate	206	—	Phenyl benzoate	—	69
Ethyl octanoate	208	—	Methyl 3-nitrobenzoate	—	78
Ethyl cyanoacetate	208	—	Methyl 4-bromobenzoate	—	81
Ethyl benzoate	212	—	Ethyl 4-aminobenzoate	—	89
Benzyl acetate	217	—	Methyl 4-nitrobenzoate	—	96
Diethyl succinate	217	—	Dimethyl fumarate	—	102
Diethyl fumarate	219	—	Cholesterol acetate	—	114
Methyl phenylacetate	220	—	Ethyl 4-hydroxybenzoate	—	116

Procedures for Preparing Derivatives

CAUTION

Some of the chemicals used in preparing derivatives are suspected carcinogens. Before beginning any of these procedures, consult the list of suspected carcinogens in Technique 1, Section 1.4 Exercise care in handling these substances.

ALDEHYDES AND KETONES

Semicarbazones

Place 0.5 mL of a 2M stock solution of semicarbazide hydrochloride (or 0.5 mL of a solution prepared by dissolving 1.11 g of semicarbazide hydrochloride [$MW = 111.5$] in 5 mL of water) in a small test tube. Add 0.15 g of the unknown compound to the test tube. If the unknown does not dissolve in the solution or if the solution becomes cloudy, add enough methanol (maximum of 2 mL) to dissolve the solid and produce a clear solution. If a solid or cloudiness remains after adding 2 mL of methanol, do not add any more methanol and continue this procedure with the solid present. Using a Pasteur pipette, add 10 drops of pyridine and heat the mixture in a hot water bath (about 60°C) for about 10–15 minutes. By that time, the product should have begun to crystallize. If the product does not crystallize, evaporate $\frac{1}{2}$ the volume of methanol. Collect the product by vacuum filtration. The product can be recrystallized from ethanol if necessary.

Semicarbazones (Alternative Method)

Dissolve 0.25 g of semicarbazide hydrochloride and 0.38 g of sodium acetate in 1.3 mL of water. Then dissolve 0.25 g of the unknown in 2.5 mL of ethanol. Mix the two solutions together in a 25-mL Erlenmeyer flask and heat the mixture to boiling for about 5 minutes. After heating the mixture, place the reaction flask in a beaker of ice and scratch the sides of the flask with a glass rod to induce crystallization of the derivative. Collect the derivative by vacuum filtration and recrystallize it from ethanol.

2,4-Dinitrophenylhydrazones

Place 2 mL of a solution of 2,4-dinitrophenylhydrazine (prepared as described for the classification test in Experiment 57D) in a test tube and add 0.15 g of the unknown compound. If the unknown is a solid, it should be dissolved in the minimum amount of 95% ethanol or 1,2-dimethoxyethane before it is added. If crystallization is not immediate, gently warm the solution for a minute in a hot water bath (90°C) and then set it aside to crystallize. Collect the product by vacuum filtration.

It is a good idea to make sure the derivative is dry. To dry a “wet” derivative, rinse the solid on the filter sequentially with a few drops of ethanol, followed by a

few drops of ether and finally with a few drops of hexane. You should avoid using too much of these solvents to avoid dissolving the derivative.

CARBOXYLIC ACIDS

Working in a hood, place 0.50 g of the acid and 2 mL of thionyl chloride into a small round-bottom flask. Add a magnetic stir bar, and attach a water-jacketed condenser and a drying tube packed with calcium chloride to the flask. While stirring, heat the reaction mixture to boiling for 30 minutes on a hot plate. Allow the mixture to cool to room temperature. Use this mixture to prepare the amide, anilide, or *p*-toluidide derivatives by one of the following three procedures.

Amides

Working in a hood, add the thionyl chloride/carboxylic acid mixture dropwise from a Pasteur pipette into a beaker containing 5 mL of ice-cold concentrated ammonium hydroxide. The reaction is very exothermic. Stir the mixture vigorously after the addition for about 5 minutes. When the reaction is complete, collect the product by vacuum filtration and recrystallize it from water or from water-ethanol, using the mixed-solvents method (see Technique 11, Section 11.10).

Anilides

Dissolve 0.5 g of aniline in 10 mL of methylene chloride in a 50-mL Erlenmeyer flask. Using a Pasteur pipette, carefully add the mixture of thionyl chloride/carboxylic acid to this solution. Warm the mixture for an additional 5 minutes on a hot plate, add a magnetic stir bar, and stir the mixture for 20 minutes at room temperature. Then transfer the methylene chloride solution to a small separatory funnel and wash it sequentially with 2.5 mL of water, 2.5 mL of 5% hydrochloric acid, 2.5 mL of 5% sodium hydroxide, and a second 2.5-mL portion of water (the methylene chloride solution should be the bottom layer). Dry the methylene chloride layer over a small amount of anhydrous sodium sulfate. Decant the methylene chloride layer away from the drying agent into a small flask and evaporate the methylene chloride on a warm hot plate in the hood. Use a stream of air or nitrogen to speed up the evaporation. Recrystallize the product from water or from ethanol-water, using the mixed-solvents method (see Technique 11, Section 11.10).

p-Toluidides

Use the same procedure as that described in preparing anilides, but substitute *p*-toluidine for aniline.

PHENOLS

α -Naphthyl-urethanes

Follow the procedure given later for preparing phenylurethanes from alcohols, but substitute α -naphthylisocyanate for phenylisocyanate.

Bromo Derivatives

First, if a stock brominating solution is not available, prepare one by dissolving 0.75 g of potassium bromide in 5 mL of water and adding 0.5 g of bromine. Dissolve 0.1 g of the phenol in 1 mL of methanol or 1,2-dimethoxyethane; then add 1 mL of water. Add 1 mL of the brominating mixture to the phenol solution and swirl the mixture vigorously. Then continue adding the brominating solution dropwise while swirling, until the color of the bromine reagent persists. Finally, add 3–5 mL of water and shake the mixture vigorously. Collect the precipitated product by vacuum filtration and wash it well with water. Recrystallize the derivative from methanol-water, using the mixed-solvents method (see Technique 11, Section 11.10).

AMINES

Acetamides

Place 0.15 g of the amine and 0.5 mL of acetic anhydride in a small Erlenmeyer flask. Heat the mixture for about 5 minutes; then add 5 mL of water and stir the solution vigorously to precipitate the product and hydrolyze the excess acetic anhydride. If the product does not crystallize, it may be necessary to scratch the walls of the flask with a glass rod. Collect the crystals by vacuum filtration and wash them with several portions of cold 5% hydrochloric acid. Recrystallize the derivative from methanol–water, using the mixed-solvents method (see Technique 11, Section 11.10).

Aromatic amines, or those amines that are not very basic, may require pyridine (2 mL) as a solvent and a catalyst for the reaction. If pyridine is used, a longer period of heating is required (up to 1 hour), and the reaction should be carried out in an apparatus equipped with a reflux condenser. After reflux, the reaction mixture must be extracted with 5–10 mL of 5% sulfuric acid to remove the pyridine.

Benzamides

Using a centrifuge tube, suspend 0.15 g of the amine in 1 mL of 10% sodium hydroxide solution and add 0.5 g of benzoyl chloride. Cap the tube and shake the mixture vigorously for about 10 minutes. After shaking the mixture, add enough dilute hydrochloric acid to bring the pH of the solution to pH 7 or 8. Collect the precipitate by vacuum filtration, wash it thoroughly with cold water, and recrystallize it from ethanol–water, using the mixed-solvents method (see Technique 11, Section 11.10).

Benzamides (Alternative Method)

In a small round-bottom flask, dissolve 0.25 g of the amine in a solution of 1.2 mL of pyridine and 2.5 mL of toluene. Add 0.25 mL of benzoyl chloride to the solution, and heat the mixture under reflux for about 30 minutes. Pour the cooled reaction mixture into 25 mL of water, and stir the mixture vigorously to hydrolyze the excess benzoyl chloride. Separate the toluene layer and wash it, first with 1.5 mL of water, and then with 1.5 mL of 5% sodium carbonate. Dry the toluene over granular anhydrous sodium sulfate, decant the toluene into a small Erlenmeyer flask, and remove the toluene by evaporation on a hot plate in the hood. Use a stream of air or nitrogen to speed up the evaporation. Recrystallize the benzamide from ethanol or ethanol–water, using the mixed-solvents method (see Technique 11, Section 11.10).

Methiodides

Mix equal-volume quantities of the amine and methyl iodide in a small round-bottom flask (about 0.25 mL is sufficient) and allow the mixture to stand for several minutes. Then heat the mixture gently under reflux for about 5 minutes. The methiodide should crystallize on cooling. If it does not, you can induce crystallization by scratching the walls of the flask with a glass rod. Collect the product by vacuum filtration and recrystallize it from ethanol or ethyl acetate.

ALCOHOLS

3,5-Dinitrobenzoates

Liquid Alcohols

Dissolve 0.25 g of 3,5-dinitrobenzoyl chloride in 0.25 mL of the alcohol and heat the mixture for about 5 minutes.¹ Allow the mixture to cool and add 1.5 mL of a

¹3,5-Dinitrobenzoyl chloride is an acid chloride and hydrolyzes readily. The purity of this reagent should be checked before its use by determining its melting point (mp 69–71 C). When the carboxylic acid is present, the melting point will be high.

5% sodium carbonate solution and 1 mL of water. Stir the mixture vigorously and crush any solid that forms. Collect the product by vacuum filtration, and wash it with cold water. Recrystallize the derivative from ethanol–water, using the mixed-solvents method (see Technique 11, Section 11.10).

Solid Alcohols

Dissolve 0.25 g of the alcohol in 1.5 mL of dry pyridine and add 0.25 g of 3,5-dinitrobenzoyl chloride. Heat the mixture under reflux for 15 minutes. Pour the cooled reaction mixture into a cold mixture of 2.5 mL of 5% sodium carbonate and 2.5 mL of water. Keep the solution cooled in an ice bath until the product crystallizes, and stir it vigorously during the entire period. Collect the product by vacuum filtration, wash it with cold water, and recrystallize it from ethanol–water, using the mixed-solvents method (see Technique 11, Section 11.10).

Phenylurethanes

Place 0.25 g of the *anhydrous* alcohol in a dry test tube and add 0.25 mL of phenylisocyanate (*α*-naphthylisocyanate for a phenol). If the compound is a phenol, add 1 drop of pyridine to catalyze the reaction. If the reaction is not spontaneous, heat the mixture in a hot water bath (90°C) for 5–10 minutes. Cool the test tube in a beaker of ice, and scratch the tube with a glass rod to induce crystallization. Decant the liquid from the solid product or, if necessary, collect the product by vacuum filtration. Dissolve the product in 2.5–3 mL of hot ligroin or hexane, and filter the mixture by gravity (preheat funnel) to remove any unwanted and insoluble diphenylurea present. Cool the filtrate to induce crystallization of the urethane. Collect the product by vacuum filtration.

ESTERS

We recommend that esters be characterized by spectroscopic methods whenever possible. A derivative of the alcohol part of an ester can be prepared with the following procedure. For other derivatives, consult a comprehensive textbook. Several are listed in Experiment 57I.

3,5-Dinitrobenzoates

Place 1.0 mL of the ester and 0.75 g of 3,5-dinitrobenzoic acid in a small round-bottom flask. Add 2 drops of concentrated sulfuric acid and a magnetic stir bar to the flask and attach a condenser. If the boiling point of the ester is above 150°C, heat at reflux while stirring for 30–45 minutes. If the boiling point of the ester is above 150°C, heat the mixture at about 150°C for 30–45 minutes. Often the mixture will turn black during the heating process. You should continue the reflux regardless. Cool the mixture, and transfer it to a small separatory funnel. Add 10 mL of ether. Extract the ether layer 2 times with 5 mL of 5% aqueous sodium carbonate (save the ether layer). Wash the organic layer with 5 mL of water, and dry the ether solution over magnesium sulfate. Evaporate the ether in a hot water bath in the hood. Use a stream of air or nitrogen to speed the evaporation. Dissolve the residue, usually an oil, in 2 mL of boiling ethanol and add water dropwise until the mixture becomes cloudy. Cool the solution to induce crystallization of the derivative.

Preparation of a Solid Carboxylic Acid from an Ester

An excellent derivative of an ester can be prepared by a basic hydrolysis of an ester when it yields a solid carboxylic acid. A procedure is provided in Experiment 57I. Melting points for solid carboxylic acids are included in the Carboxylic Acids Table in Appendix 1.

Index of Spectra

Infrared Spectra

Anisole	904
Benzaldehyde	310, 907
Benzamide	911
Benzil	312
Benzilic acid	315
Benzocaine	408
Benzoic acid	326, 909
Benzoin	310
Benzonitrile	906
Borneol	297
1-Bromobutane	202
<i>n</i> -Butyl bromide	202
<i>n</i> -Butylamine	905
Caffeine	106
Camphor	297, 909
Carbon disulfide	892
Carbon tetrachloride	889
Carvone	133
Chloroform	890
Decane	899
1,2-Dichlorobenzene	902
<i>N,N</i> -Diethyl- <i>m</i> -toluamide	421
6-Ethoxycarbonyl-3,5-diphenyl-2-cyclohexenone	364
Isoborneol	298
Isopentyl Acetate	114, 910
Limonene	133
Mesityl oxide	908
2-Methyl-4-nitro-1-(pent-1-yn-1-yl)benzene	339
Methyl benzoate	910
Methyl isopropyl ketone	895
Methyl <i>m</i> -nitrobenzoate	385
4-Methylcyclohexanol	212, 904
4-Methylcyclohexene	211, 900
Mineral oil	888
2-Naphthol	903
Nitrobenzene	906
Nonanal	907
<i>cis</i> -Norborene-5,6- <i>endo</i> -dicarboxylic anhydride	458, 911
Nujol	888
<i>t</i> -Pentyl chloride	204
Polystyrene	893

Styrene	901
Sulfanilamide	430
Triphenylmethanol	324

¹H NMR Spectra

Benzocaine	408
Benzyl acetate (60-MHz spectrum)	922
Benzyl acetate (300-MHz spectrum)	921
Borneol	299
Camphor	298
(-)-Carvone	134
<i>N,N</i> -Diethyl- <i>m</i> -toluamide	422
6-Ethoxycarbonyl-3,5-diphenyl-2-cyclohexenone	365
Ethyl 4-(pent-1-yn-1-yl)benzoate	338
Ethyl 3-hydroxybutanoate	269
Eugenol	123
1-(Hex-1-yn-1-yl)-4-nitrobenzene	941
1-(4-(Hex-1-yn-1-yl)phenyl)ethanone	337
Hexanol (1-)	338
(<i>E</i>)-4-(4-Hydroxy-3-methoxyphenyl)-2-buten-1-ol	347
9-(Hydroxymethyl)-13-methyl-10,11-dihydro-9 <i>H</i> -9,10-[3,4]epipyrrroloanthracene-12,14(13 <i>H</i> ,15 <i>H</i>)dione	944
Isoborneol	299
(+)-Limonene	134
Naproxen	337
Phenylacetone	561
1,1,2-Trichloroethane	915
Vegetable oil	927
Vinyl acetate	931

¹³C NMR Spectra

Borneol	300
Camphor	300
Carvone	135
Cyclohexene	959
Cyclohexanol	960
Cyclohexanone	959
2,2-Dimethylbutane	958
Ethyl phenylacetate	956
Isoborneol	301
1-Propanol	957
Toluene	961

Mass Spectra

Acetophenone	982
Benzaldehyde	981
Bromoethane	974
1-Bromohexane	985
Butane	976

1-Butanol	980
2-Butanone	982
1-Butene	979
Chloroethane	973
Cyclopentane	978
Dopamine	970
Methyl butanoate	984
Propanoic acid	983
Toluene	979
2,2,4-Trimethylpentane	977

Ultraviolet-Visible Spectra

Benzophenone	466
Naphthalene	466

Mixtures

Borneol and isoborneol	302
1-Bromobutane and 1-Chlorobutane	196
<i>tert</i> -Butyl bromide and <i>tert</i> -Butyl chloride	196

Index

A

- Ab initio* calculations, 172
- Acetaminophen, 77
 - crude acetaminophen isolation, 83, 84–85
 - crystallization, 83–84, 85
 - decolorization, 83, 85
 - heating, 82, 84
 - microscale procedure, 82–84
 - reaction mixture, 82, 84
 - semimicroscale procedure, 84–86
 - yield calculation, 84, 85–86
- Acetylsalicylic acid, 69–72, 77
 - aspirin tablets, 72
 - crystallization of, 71
 - ferric chloride test, 71
 - melting point, 71
 - procedure, 70–72
 - vacuum filtration, 71
- Achiral environment, 270
- Acid chloride, 419–420
- Acids, 908
- Acid-washed alumina, 809
- Activators, 377
- Acylglycerols, 213
- Addition funnel, 654
- Addition polymers, 433–434, 435t
- Addition reaction, 501
- Air peak, 855
- Alarm pheromones, 410, 413
- Alcohol, 518–523, 903
 - cerium (IV) test, 519–521
 - chromic acid test, 521–522
 - classification tests, 519–523
 - derivatives, 523
 - iodoform test, 522
 - Lucas test, 521
 - spectroscopy, 523
- Aldehyde, 504–510, 905
 - chromic acid test, 507–508
 - classification tests, 505–510
 - derivatives, 510
 - 2,4-dinitrophenylhydrazine, 505–506
 - ferric chloride test, 509–510
 - iodoform test, 508–509
 - spectroscopy, 509
 - Tollens test, 506–507
- Aldehyde condensation reaction, 582–585
- Aldehyde enigma, 566–568
 - aqueous layer, 567
 - organic layer, 567
 - procedure, 566–567
- Aldol condensation reaction, 349–353, 569, 582–585
 - green enantioselective, 354–360
 - molecular modeling, 351–352
 - procedure, 350–351
- Alembic apparatus, 756, 757f
- Alizarin, 386
- Alkaloids, 94
- Alkanes, 899
- Alkenes, 899
- Alkyl halides, 184–188
 - molecular modeling, 187–188
 - procedure, 186–187
- Alkynes, 901
- Allergens, 545
- Allinger, Norman, 161
- Aluminum block, 2–3, 3f
- Ambident nucleophile, 352
- Amides, 418, 910
 - isolation of, 420–421
 - preparation of, 420
- Amines, 515–518, 903–904
 - Hinsberg test, 517
 - nitrous acid test, 516
- Amoore, J.E., 125
- Analgesic drug
 - active ingredient isolation, 78–80
 - thin-layer chromatography, 89–93
 - commercial analgesics analysis, 93
 - development chamber, 92
 - first reference plate, 91–92
 - initial preparation, 91–92

- Analgesic drug (*Continued*)
 iodine analysis, 93
 plate development, 92
 UV visualization, 92–93
- Analgesic drug, active ingredient isolation, 77–80
- Analgesics
 caffeine and, 75
 column chromatography, 79
 history of, 73
 solvent evaporation, 79–80
 vacuum filtration, 80
- Analytical balance, 859
- Anhydrides, 911
- Anthracene-9-methanol, 459–461
- Antihistamine drugs analysis, 545–547
- APC tablet, 74
- Aromatic compounds
 nitration, 274–277
 procedure, 379–380
 relative reactivities of, 377–380
- Aromaticity tests, 503–504
- Aromatic rings, 900–901
- Artificial pineapple flavor, 109t
- Aspirin. *See also* Acetylsalicylic acid
 chemical structure of, 67
 description of, 66
 history of, 66
- Atom economy, 260
- Atomic orbital, linear combination, 169
- ATR. *See* Attenuated total reflectance
- Attenuated total reflectance (ATR), 885
- Automatic pipettes, 6–7, 6f
- Axel, Richard, 126
- Azeotropes
 distillation application, 780–782, 781f, 782f
 generalizations, 779
 maximum boiling point, 779t
 minimum boiling point, 777–779, 777f, 778t
 nature of, 776f
- Azeotropic distillation, 780–782, 781f, 782f
- B**
- Back extraction, 735
- Base peak, 970
- Basic alumina, 808
- Basis-set orbital
 Gaussian-type orbital, 170
 polarization basis sets, 171
 Slater-type orbital, 170
 split-valence basis sets, 170–171
- Benzaldehyde conversion, 304–316
- Benzilic acid preparation, 313–316
 procedure, 314–315
- Benzil preparation
 procedure, 312–313
- Benzocaine, 406–409
- Benzoic acid, 324–327
 crystallization, 326
 dry ice addition, 324–325
 hydrolysis, 325
 product isolation, 325
 spectroscopy, 326
- Benzoin preparation
 procedure, 309
 by thiamine catalysts, 305–311
- Benzophenone, photoreduction of
 energy transfer, 466–467
 fluorescence, 463
 internal conversion, 466
 intersystem crossing, 463
 phosphorescence, 464
 procedure, 468
 radiationless transition, 464
 singlet state, 463
- β -Benzopinacolone, synthesis of, 469–470
- Beriberi, 305
- Binding energy, 172
- Biodiesel, 250–251
 analysis of, 255
 calorimetry, 255
 from coconut oil, 254
 from other oils, 255
 spectroscopy, 255
- Biofuels
 biodiesel, 250–251
 ethanol, 248–249
 nature of, 248
- Bioluminescence, 471
- Biphenyl, 319
- Boiling point
 determination of
 difficulties for, 750
 microscale inverted capillary method, 748–750
 semimicroscale direct method, 746–747, 747f
 semimicroscale inverted capillary method, 748
 digital thermometers, 752–753, 753f
 glass thermometers, 750–752, 752f
 nature of, 745, 745f, 746f
 stem corrections, 750–752, 752f
- Boiling-point determination, 62–65
 analysis, 64
 report, 64
- Boiling stone, 652–653
- Bond-density surface, 174
- Borneol, oxidation of, 289–290
 chromium removal, 294
 reaction setup, 293–294
 solvent removal, 294
 sublimation procedure, 295–296
 waste disposal, 294–295
- Büchner funnel, 673
- Buck, Linda, 126f
- C**
- C-4 or C-5 acetate ester, 528–531
- Caffeine
 carbon dioxide decaffeination process, 96
 content per liquid, 95, 96
 decaffeinated coffee, 96
 direct contact decaffeination, 96
 as diuretic, 95
 effect of, 94–95
 espresso coffee, 95
 history of, 94

- isolation from coffee, 98–106
- isolation from tea, 98–106
 - drying, 101–102
 - evaporation, 102
 - extraction, 101–102
 - sublimation, 102–103
 - tea solution, 101
- solid extraction, 103–106
- as stimulant, 95
- in tea, 95
- water process, 96
- Camphor, reduction of, 290–292
 - analysis, 296–297
 - gas chromatography, 297
 - isolation, 296–297
 - molecular modeling, 302–303
 - reduction, 296
- Caraway oil, 129–136
 - carvones analysis
 - boiling point, 133
 - gas chromatography, 132, 135–136
 - infrared spectroscopy, 132
 - nuclear magnetic resonance spectroscopy, 132–133
 - odor, 132
 - polarimetry, 132
 - refractive index, 133
 - structure of, 130
- Carbon-13 nuclear magnetic resonance spectroscopy, 952–968
 - aromatic ring compounds, 960–962, 961f
 - chemical shifts, 953–954, 953f
 - proton decoupling, 956
 - sample preparation, 952–953
 - sample spectra, 958–959, 958f
- Carbonyl compounds, 905
- Carboxylic acids, 510–512
- Carotenoids, 139
- Catalogs, 625–631
- Cellulosic ethanol, 249
- Chalcones
 - cyclopropanation, 578–581
 - green epoxidation, 574–577
 - substituted chalcones synthesis, 569–573
 - crude product isolation, 572
 - crystallization, 572
 - laboratory report, 572
 - procedure, 571–573
- Chemical literature guide, 987–1000
 - advanced laboratory techniques, 992
 - advanced textbooks, 989–990
 - Beilstein* and *Chemical Abstracts*, 994–996
 - computer online searching, 996–997
 - current interest, 998
 - general synthetic methods, 988
 - literature search guide, 999
 - organic qualitative analysis, 993–994
 - physical constants, 987–988
 - reaction mechanisms, 992–993
 - researching literature, 988–989
 - scientific journals, 997–998
 - specific synthetic methods, 990–992
 - spectra collections, 989
- Chemical reagent, 306
- Chemical shift, 919–920, 924t
- Chemical shift reagents, 940–941
- Chemiluminescence, 471
- Chemotherapy, 423
- Chiral environment, 270
- Chiral recognition, 129
- Chlorophyll
 - chlorophyll *a*, 143
 - chlorophyll *b*, 143
 - definition of, 142
 - pheophytin *a*, 143
 - pheophytin *b*, 143
- Chlorophyll, definition of, 142
- Chloroplasts, 142
- Chromatography, 45–53, 88
- Clove essential oils
 - microscale procedure
 - drying, 122
 - evaporation, 122
 - extraction, 122
 - steam distillation, 121–122
 - yield determination, 122
 - nature of, 120
 - semimicroscale procedure, 123–124
 - steam distillation, 120–124
- Cocaine, 402
- Coffee
 - caffeine isolation, 101–103
 - solid extraction, 103–106
 - analysis, 105–106
 - caffeine isolation, 105
 - caffeine removal, 105
 - solution filtering, 104–105
 - solution preparation, 104
 - sublimation, 105
- Cold pressing, 214
- Colorants, 399
- Column chromatography, 50–53, 51f
 - adsorbents, 808–809
 - column packing, 815–816
 - microscale methods, 816–817, 817–820
 - monitoring, 822–823
 - problems with, 815–816
 - sample application, 820
 - semimicroscale methods, 817–820
 - separated compounds, 823
 - tailing, 823, 824f
 - column running, 51–52
 - decolorization by, 824
 - definition of, 808
 - elution techniques, 821
 - flash chromatography, 826, 826f
 - gel chromatography, 824–825, 825f
 - interactions, 809–811, 809f, 810f, 811f
 - preparation, 51
 - reservoirs, 821–822, 822f
 - separation parameters
 - adsorbents, 812–813, 812t
 - column size, 815, 815t
 - flow rate, 815
 - solvents, 813–814, 813t, 814t
 - types of, 808

Computational chemistry, 158, 176–181
 terms for, 168–169
 Condensation polymers, 434–439, 436t
 Conformation, 160
 Conical reaction vials, 3–5, 4f
 Correlation chart, 896, 897t
 Craig tube, 20, 83–84, 675, 675f, 676f, 705f
 Critical thinking application, 28–31, 39–41
 Crystallization, 20–31, 69
 common solvents, 711t
 decolorization, 712–713
 inducing of, 713
 microscale crystallization, 24–27, 705f, 708t
 Craig tube, 705f
 decolorization, 712
 drying the crystals, 707
 inducing of, 713
 insoluble impurities removal, 706–707
 solid crystallization, 707
 solid dissolving, 706
 mixed solvents, 715
 semimicroscale crystallization, 21–24, 699–704,
 700f, 708t
 crystal collection, 704
 decolorization, 712
 dissolving, 701
 inducing of, 713
 insoluble impurities and, 701–704, 703f
 solid crystallization, 704
 solids purification, 696–717
 solubility, 696–697, 697f
 theory of, 698–699
 solvent selection, 708–710
 solvent testing, 710–711
 Current implementations, 161–162
 Cyano group detection, 497–498
 Cyclopentadiene, 455–458

D

Daruma, 94
 Density, 754
 Density-electrostatic potential, 174
 carbocations, 180
 carbonyl reactivities, 180–181
 carboxylic acids, acidities of, 179–180
 LUMO maps, 180–181
 maps, 179–180
 Density surface, mapping properties of, 174–175
 Detergent, 399
 Dewar flask, 792
 Diastereotopic, 461, 548
 Diels–Alder reaction, 449, 455–458, 459–461
 Dienophile, 449
 Dimethylcyclohexanes, 166
 1,4-Diphenyl-1, 3-butadiene, 369–376
 preparation of, 372–374
 potassium phosphate, 374–376
 sodium ethoxide, 372–374
 Disparlure, 412
 Dispensing pumps, 7–8, 7f

Distribution coefficient, 719–721
 Domagk, Gerhard, 423
 Drug identification, 87–88
 Drying agents, 718–744
 Duisberg, Carl, 74
 Dyes, 399–401

E

Ebulliator tube, 788, 788f
 Eckert, Charles, 261
 Ehrlich, Paul, 423
 Electron-density surface, 174
 Electrophilic addition, 381
 Electrophilic aromatic substitution reactions, 382
 Electrophilic reagents, 381
 Elements, test for, 494–500
 Eluates, 810
 Eluents, 810
 Elutants, 810
 Enantioselective process, 264
 Enantiospecific process, 264
 Enol dianion, 473
 Enzyme, competitive inhibitor of, 425
 Equal, 447
 Essential oils
 spectroscopy, 123–124
 steam distillation of, 120–124
 Esters, 523–526, 909
 basic hydrolysis, 524–525
 definition of, 107
 derivatives, 526
 ferric hydroxamate test, 524–525
 flavors, 107
 spectroscopy, 525–526
 Ethanol, 149–151
 distillate analysis, 156
 fermentation, 154
 fermentation apparatus, 154
 fractional distillation, 155
 from sucrose, 152–156
 Ethers, 903
 Ethyl acetoacetate
 chiral reduction of, 264–273
 alcohol isolation, 266–267
 gas chromatography, 268
 infrared spectroscopy, 267–268
 polarimetry, 268
 yeast reduction, 266
 Ethyl (S)-3-hydroxybutanoate, 269–273, 272f
 Extraction, 9–10, 32–41
 caffeine extraction, 33–35
 distribution coefficient, 719–721, 720f, 721f
 drying agents, 728–732, 729t
 emulsions, 732–733
 flowchart usage, 735–738, 737f
 general process of, 718, 719f
 liquid extraction, 738–739, 739f
 method selection, 721–722, 722t
 microscale extraction
 conical vial, 722–726, 723f, 724f, 725f

- drying procedure, 730–732, 731t
 - separatory funnel, 726–727, 726f, 727f
 - neutral compound, 37–39
 - organic layer, 727–728
 - purification methods, 733–735, 734f
 - separation methods, 733–735, 734f
 - solid phase extraction, 739–742, 740f, 741f
 - Extractions, 718–744
- F**
- Fats, 213–217
 - common fatty acids, 214t
 - fat replacement, 216
 - fatty acid composition, 215t
 - rendering of, 214
 - trans*-fatty acids, 216
 - FDA. *See* Federal Drug Administration; Food and Drug Administration
 - Federal Drug Administration (FDA), 437
 - Fermentation chemistry, 149–151
 - Ferrous hydroxide test, 497
 - Filtration, 667–677
 - aspirator, 674–675, 674f
 - centrifugation, 676
 - Craig tube, 675, 675f, 676f
 - filtering media, 673
 - filter paper, 671–672, 671t
 - filter-tip pipette, 675, 675f
 - gravity filtration
 - decantation, 671
 - filter cones, 668, 668f
 - filtering pipettes, 670–671, 670f
 - fluted filters, 668–670, 669f
 - methods, 667t
 - vacuum filtration, 672–673, 672f
 - Fireflies, 471–473
 - Flame front, 236
 - Food and Drug Administration (FDA), 73, 76, 95
 - Force constant, 160
 - Four-cycle engine, 235
 - Fractional distillation, 54–61, 758
 - apparatus, 56–57, 57f
 - apparatus for, 769f, 776f
 - azeotropes, 768–784
 - column efficiency, 773
 - distillation, 57–59, 58f
 - distillation curve, 59
 - gas chromatography, 59–60
 - instructor notes, 55–56
 - methods for, 775–776
 - microscale procedure, 60–61
 - practice for, 775–776
 - Raoult's Law, 771–772, 772f
 - types of columns, 773–775, 774f, 775f
 - vapor-liquid composition diagrams, 770–771, 771f
 - vs.* simple distillation, 768–770
 - Fragment ion peak, 970
 - Friedel-Crafts acylation, 537–544
 - carbon-13 NMR spectrum, 538–539
 - infrared spectrum, 538
 - instructions, 539
 - procedure, 541–543, 541f
 - proton-NMR spectrum, 538
 - report, 543
 - Frozen core approximation, 171–172
- G**
- Gas chromatography, 847–866
 - advantages of, 853
 - apparatus, 847, 848f
 - chiral stationary phases, 855–857, 856f, 857f
 - columns, 849–851
 - capillary columns, 850–851
 - monitoring of, 853–855, 854f, 855f
 - packed columns, 849–850
 - data tables, 862f, 863f
 - component percentage, 864
 - response factor application, 862–863
 - GC-MS, 864–865
 - liquid phases, 849t
 - moving gas phase, 847
 - qualitative analysis, 857
 - quantitative analysis, 859–861, 859f, 860f
 - retention time, 855
 - sample collection, 858–859, 858f
 - separation factors, 852–853
 - separation principles, 851–852, 852f
 - solid supports, 849t
 - stationary liquid phase, 847
 - tailing, 857, 857f
 - Gas chromatography-mass spectrometry (GC-MS), 545, 864–865
 - Gasolines
 - gas-chromatographic analysis of, 243–247
 - major components, 246
 - procedure
 - analysis, 245
 - gasoline samples, 245
 - oxygenated fuel reference mixture, 245
 - reference mixture, 245
 - GC-MS. *See* Gas chromatography-mass spectrometry
 - General unknown, 480
 - Global minimum, 160
 - Glycerides, 213
 - Gossypure, 412
 - Gradient elution systems, 845
 - Graduated pipettes, 8
 - Grain alcohol (ethyl alcohol), 149
 - Graphic models, 173–174
 - Green chemistry, 258–263
 - atom economy, 260
 - biomass, 259
 - ionic liquids, 261
 - principles of, 259
 - Green enantioselective aldol condensation reaction, 354–360
 - Greenhouse effect, 241
 - Grignard reagent, 320–321, 322
 - Grubbs-catalyzed metathesis, 338–344
 - Gyplure, 412

- H**
- Halide, tests for, 494–497, 499
- Halides, 912
- Handbooks, 625–631
Aldrich handbook of fine chemicals, 630–649
CRC handbook, 625–627, 626t, 627t
Lange's handbook of chemistry, 627–629, 628t, 649t
Merck index, 629–648
- Hard pesticides, 450
- Hartrees, 173
- Heating/cooling methods, 640–646, 640f, 641f
 aluminum block, 641f
 cold baths, 645
 dial thermometers, 642f
 flames, 644–645
 sand bath, 643, 643f
 steam baths, 645–646, 645f
 water bath, 644, 644f
- Heating methods, 2–3
 aluminum block, 2–3, 3f
 sand baths, 3
 water bath, 3
- Heat of formation, 172–173
- Heats of reaction, 178–179
- Helianthin, 392
- Henna, 386
- HETP, 775
- Hickman head, 759, 759f, 781
- High-performance liquid chromatography (HPLC), 842–846, 843f
 adsorbents, 843–844
 C-18 column, 844
 columns, 843–844
 dimensions, 844–845
 data presentation, 846
 detectors, 846
 preparative column, 845
 reversed-phase chromatography, 844
 solvents, 845
- Hinsberg test, 517
- Hirsch funnel, 23, 85, 673, 699
- Hofmann, Felix, 66
- Hot pressing, 214
- HPLC. *See* High-performance liquid chromatography
- Hubbert, Marion King, 240
- Humectant, 399
- Hund's rule, 463
- Hydrolysis, 69
- Hydrophobic effect, 459
- I**
- Ibuprofen, 75, 91
- Immiscible solvents, 35–36
- Indigo, preparation of, 396–398
 reduction of, 398
 synthesis of, 397
 UV characterization, 397–398
- Indoxyl, 386
- Inert compounds, 493
- Infrared spectroscopy, 62–65, 880–913
 acids, 908
 alcohol, 903
 aldehydes, 905
 alkanes, 899
 alkenes, 899
 alkynes, 901
 amides, 910
 amines, 903–904
 analysis, 64
 anhydrides, 911
 aromatic rings, 900–901
 calibration, 892–893
 carbonyl compounds, 905
 correlation chart, 896, 897t
 esters, 909
 ethers, 903
 halides, 912
 infrared spectrum, 63–64
 ketones, 907
 nitro compounds, 905
 phenols, 903
 report, 64
 sample preparation
 introduction, 881–882
 liquid samples, 882–883, 883f, 884–885, 884f
 solid samples, 885–893, 889f, 890f
 spectrum recording, 892
 uses of, 893
 vibration modes, 894
- Injection port, 851
- Insecticides, 449–453
 alternatives to, 452–453
- Insect repellent, 414–415
- Invertase, 152, 152
- Invert sugar, 152
- Iodosubstituted aromatic compounds
 microwave technology, 335–336
 procedure
 column chromatography, 331–332
 product analysis, 332–335, 333f, 334f, 335f
 product isolation, 332
 reaction mixture, 331
 Sonogashira coupling of, 328–344
- Ion-exchange chromatography, 844
- Ionic liquids, 261
- Ionization chamber, 969
- Isomerism, 177–178
- Isopentyl acetate (banana oil)
 microscale procedure
 apparatus, 111
 distillation, 112
 preparation, 111–112
 reflux, 112
 workup, 112
 yield determination, 112
 semimicroscale procedure
 apparatus, 113
 distillation, 114
 drying, 113
 extractions, 113

reaction mixture, 113
 reflux, 113
 yield determination, 114
 Isoprene rule, 117
 Isoprene units, 117

J

Juvenile hormone, 452

K

KBr pellets, 886–888, 886f
 Ketones, 504–510, 907
 chromic acid test, 507
 classification tests, 505–510
 derivatives, 510
 ferric chloride test, 509
 iodoform test, 508
 spectroscopy, 509–510
 Tollens test, 506
 Ketoprofen, 76
 Kuo P'ou, 94

L

Laboratory glassware, 617–624
 assembling, 620
 cleaning of, 617–618
 drying of, 618
 equipment description, 621–642, 621f, 640f, 641f, 642f
 etching glassware, 620
 ground-glass joints, 618–619, 619f
 Laboratory notebook, 609–616
 actual yield, 615
 advance preparation, 610
 calculations, 612
 format, 610–616
 laboratory records, 612
 laboratory reports, 615
 limiting reagent, 612
 percentage yield, 615
 preparative experiments, 610
 sample of, 613f, 614f
 sample submission, 615–616
 separation scheme, 611, 611f
 for side reactions, 609
 theoretical yield, 612
 weight percentage recovery, 615
 Laboratory safety, 590–608
 bottle labels, 604
 carcinogenic substances, 606–607
 clothing, 594–595
 common solvents, 604–606
 eye safety, 590
 fires, 590–591
 first aid, 595
 food, 594
 inadvertently mixed chemicals, 594
 material safety sheets, 595–602
 organic solvents, 591–592
 right-to-know laws, 595–604
 safety data sheets, 602–604
 safety guidelines, 590–595

unauthorized experiments, 594
 use of flames, 593–594
 waste disposal, 592–593

Lachrymators, 185, 346
 Leaded gasoline, 237
 Lewis acid catalyst, 378
 Limiting reactants, 632
 Liotta, Charles, 261
 Liquid, measurement of, 6–8, 6f
 dispensing pumps, 7–8, 7f
 graduated pipettes, 8
 Local anesthetics, 404f
 nature of, 402
 synthetic substitutes, 403
 Local minimum, 160
 Lucas test, 521
 Luciferase, 471
 Luciferin, 471
 Lucretius, 125, 127
 Luer lock adapter, 557
 Luminol, 474–478

M

Macroscale, 2
 Magnetic stirrers, 651
 Maleic anhydride, 455–458
 Maltase, 152
 Maltose, 152
 Mass analyzer, 969
 Mass spectrometry, 969–986
 fragmentation patterns, 972–976, 972t
 halogen detection, 972
 interpreted mass spectra, 976–984, 977f, 978f, 979f, 981f, 982f, 983f, 984f
 mass spectrum, 969–970
 molecular formula determination, 971–972
 rearrangement reactions, 985
 Mass spectrum, 969
 Mauve, 387
 McClintock, Martha, 412
 McClintock effect, 412
 McLafferty rearrangement, 985
 Measurement
 of liquid, 6–8
 of solids, 5, 5f
 Melting point, 678–686
 composition curve, 679f
 decomposition, 683–685
 discoloration, 683–685
 electrical instruments for, 682–683, 682f
 melting-point standards, 685t
 melting-point tube, 681
 mixture melting points, 680
 physical properties, 678
 shrinkage, 683–685
 softening, 683–685
 sublimation, 683–685
 theory, 679–680, 679f
 thermometer calibration, 684–685, 685f
 Thiele tube, 681, 681f
 Merck index, 16, 17, 28, 64

- Mercury thermometer, 3, 4f
- Methyl benzoate, 381–385
 molecular modeling, 384–385
 procedure, 383–384
- 4-Methylcyclohexene, 207–212
 microscale procedure, 209–210
 semimicroscale procedure, 210–212
 unsaturation tests, 212
- Methyl orange, preparation of, 392–396
 procedure, 394–395
 special instructions, 394
- Methyl stearate
 from methyl oleate, 218–222
 procedure
 apparatus, 219–220, 220f
 catalyst, 221–222
 reaction preparation, 221
- Mevalonic acid, 117
- Michael condensation reaction, 582–585
- Microscale, 2
- Microscale crystallization (Craig tube)
 crystal isolation, 26
 crystallization, 25–26
 pre-lab calculations, 24–25
 preparations, 25
 sulfanilamide, 25
- Microscale laboratory, 2–11
- Microspatulas, 10–11, 11f
- Minimization, 160
- Mixture melting point, 21, 27–28
- MM. *See* Molecular mechanics
- Molecular ion, 970
- Molecular mechanics (MM), 168
 history of, 158
 limitations of, 161
 molecular force field factors, 159t
- Molecular modeling, 158
ab initio calculation, 169
 boat conformations, 165
cis- and *trans*-2-butene, 166–167
 cyclohexane chair, 165
 methods, 168–169
n-butane conformations, 164–165
 semiempirical calculations, 169
 substituted cyclohexane rings, 166
- Molecular rotation, 869
- Molecular sieve chromatography, 825
- Moncrieff, R.W., 125
- Mordants, 386
- Mother liquor, 20, 698
- Mydriasis, 403
- N**
- Naproxen, 76
- Naproxen synthesis, 552–565
 laboratory report, 565
 procedure, 556–565
- Nucleophilic substitution, 189–197
 analysis, 194–195
 with 1-butanol or 2-butanol, 191–193
 apparatus, 191–192, 192f
 purification, 192–193
 reflux, 192
 with 2-methyl-2-propanol, 193–194
n-butyl bromide, 200–202, 201f, 202
 report, 195–196
t-pentyl chloride
 macroscale procedure, 205
 microscale procedure, 203–204
 semimicroscale procedure, 204
- Nicotine, 453
- Nitro compounds, 905
- Nitrogen test, 499
- Nitro groups, detection of, 497
- Normal phase chromatography, 103, 844
- Nuclear magnetic resonance spectroscopy, 914–951, 914f
 anisotropy, 925–926, 926f
 aromatic compounds
 monosubstituted rings, 933–935
para-disubstituted rings, 936–937
 chemical equivalence, 920–922, 920f, 921f
 chemical shift, 919–920, 924t
 chemical shift reagents, 940–941
 coupling constant, 928–930, 929t
 diastereotopic protons, 942–945
 higher field strength, 932–933
 local diamagnetic shielding, 925, 925f
 magnetic equivalence, 930–932, 931f
 sample preparation
 NMR tube, 917, 917f
 reference substance, 918–919
 spin-spin splitting, 926–928, 928f
- Nucleophilic substitution, 532–536
- Nujol mull, 888–889, 888f
- Nylon. *See* Polyamide
- O**
- Odor, stereochemical theory of, 125–128
 history of, 125
 Nobel prize theory of, 128
 nose odor receptors, 127f
 odor receptor sites, 126f
- Oils, 213–217
 fatty acid composition, 215t
- Opsin, 137
- Optically active substance, 868f
- Organozinc reactions, 338–344
 activated zinc, 347
 preparation, 347–348
 reaction, 347–348
- Organozinc reagents, 548–551
- P**
- Paint
 formulation of, 399–401
 thinner, 399
- PAN. *See* Peroxyacetyl nitrate
- Paper factor, 453
- Parameterization, 160
- Pasteur, Louis, 149, 152
- Pasteur pipette, 34, 50, 52, 101, 144

- PCC. *See* Pyridinium chlorochromate
- Penicillin, 424
- Peroxyacetyl nitrate (PAN), 239
- Petroleum/fossil fuel
- alkylation, 237
 - crude oil
 - definition, 234
 - distillation, 234t
 - crude oil definition, 234
 - fuel classification, 236
 - organic compound, octane ratings for, 236t
 - paraffins, 234
 - reforming, 237
 - straight-run gasoline, 234
- Phenols, 69, 512–515, 903
- bromine water, 514
 - cerium (IV) test, 513–514
 - classification tests, 513–515
 - ferric chloride, 513
- Phenylpropanoids
- definition of, 118–119
- Pheromones, 109, 410–415
- insect repellent, 414–415
 - other pheromones, 413
 - sex attractants, 411–413
 - types of, 410
- Phosphorescence, 464
- Photochemistry, 471–473
- Pipette calibration, 9
- Plasticizers, 432–433
- Plastics, 438t
- Polarimetry, 867–874
- optical purity, 873–874
 - polarimeter, 869–870, 869f
 - modern digital polarimeter, 872–873, 873f
 - operation of, 870–872, 871f
 - polarized light, nature of, 867–868, 867f, 868f
- Polarimetry cells, 870, 869f
- Polyamide, 443–444
- Polyesters, 442–443
- Polymers
- addition polymers, 433–434, 435t
 - chemical structure of, 431
 - condensation polymers, 434–439, 436t
 - infrared spectra of, 446–447
 - preparation of, 441–451
 - properties of, 441–451
 - thermal classification of, 432–433
 - thermoplastics, 432–433
 - types of
 - addition polymers, 431–432
 - condensation polymers, 432
 - cross-linked polymers, 432
- Polystyrene, 445–446
- Precise atomic masses, 971, 971t
- Pressure-equalizing addition funnel, 654
- Primary aromatic, 516
- Primer pheromones, 410
- Product development control, 291
- Prontosil, 423
- Propylure, 412
- Prostaglandins, 67
- Proton-coupled spectra, 955
- Proton decoupling, 956
- Pure substance, 756–757
- Purification scheme, 42–44
- Pyridinium chlorochromate (PCC), 289
- Q**
- Quadrupole broadening, 939
- Quantum mechanics, 158, 168
- Queen substance, 413
- R**
- Raoult's Law, 770, 771–772, 772f
- Reaction methods, 647–666
- adjustable metal clamps, 647, 648f
 - apparatus securing of, 647–648, 647f, 649f
 - boiling stones, 652–653
 - drying tubes, 654, 654f
 - gaseous product collection, 660–661, 660f
 - heating, 649–651, 650f
 - inert atmosphere, 655–657, 656f
 - liquid reagents, 653–654, 653f
 - microwave-assisted organic chemistry, 664–665
 - noxious gases, 657–660
 - aspirator, 659–660, 659f
 - drying tube method, 657
 - external gas traps, 658–659, 658f, 659f
 - rotary evaporator, 663, 664, 664f
 - solvent evaporation, 661–663, 662f, 663f
 - stirring method, 651–652
- Recognition pheromones, 410, 413
- Recruiting pheromones, 410
- Recrystallization, 699
- Reference cell, 891
- Reflux ring, 650, 651f
- Reformulated gasoline (RFG), 238
- Refractometry
- Abbé refractometer, 876–878, 876f, 877f
 - cleaning for, 878
 - digital refractometer, 878–879, 879f
 - refractive index, 875–876, 875f
- Regioselectivity, 177–178
- Releaser pheromones, 410
- Response factor, 860
- Retention time, 855
- Retort apparatus, 756, 757f
- Reverse phase chromatography, 103
- Reye's syndrome, 68
- RFG. *See* Reformulated gasoline
- Rhodopsin, 137, 140
- Robiquet, Jean, 94
- S**
- Salicylamide, 74, 91
- Salting out, 567
- SAM. *See* Self-assembled monolayer
- Sample cell, 891
- Schrödinger equation, solving of, 169
- Semiempirical methods, 171–172

- Semimicroscale crystallization
 components of, 23–24
 crystallization, 23
 filtration, 23
 pre-lab calculations, 22
 preparations, 22
 procedure, 22–24
 sulfanilamide, 22
- Separations, 718–744
- Separation scheme, 42–44
- Sex attractants, 411–413
- Sex pheromones, 411
- Shen Nung, 94
- Simple distillation, 54–61
 apparatus, 56–57, 57f
 digital thermometer, 764, 766f
 distillation, 57–59
 distillation curve, 59
 equipment for, 764, 765f
 evolution of, 756
 fractional distillation, 758
 gas chromatography, 59–60
 instructor notes, 55–56
 microscale equipment, 759–763, 759f, 760f, 761f, 762f
 boiling stones, 762
 condenser, 760
 degree of heating, 763
 fractions, 759–760, 761f
 Hickman head, 759, 759f
 sealed systems, 760
 temperature monitoring, 760–761, 762f, 763f
 microscale procedure, 60–61
 semimicroscale procedure, 56–60
 theory of, 756–759, 757f, 758f
 vs. fractional distillation, 768–770
- Singlet energy transfer, 466
- Single-point calculations, 171
- Single-point energy, 172
- Size-exclusion chromatography, 825
- Size-exclusion column, 844
- Smog, 239
- Sodium borohydride, 290–292
- Sodium fusion test, 498–500
- Solid phase extraction (SPE), 100, 262
- Solids
 measurement of, 5, 5f
 purification of, 696–717
 sublimation behavior of, 798
- Solubility, 12–21, 687–695
 of alcohols, 14
 of bases, 15
 behavior, prediction of, 688–692
 critical thinking application, 15–16
 definition of, 687
 immiscible pairs, 14
 miscible pairs, 14
 of organic acids, 15
 organic solvents, 692–693
 polarity, 688–691, 690t, 691f
 of solid compounds, 13
- Solubility tests, 487–493
 procedure
 in concentrated sulfuric acid, 493
 in 5% HCl, 491
 inert compounds, 493
 test compounds, 489–490
 in water, 490
- Solution cell, 890, 891f
- Solvent extraction, 214
- Solvent selection, 26–27
- SPE. *See* Solid phase extraction
- Spearmint oil, 129–136
 carvones analysis
 boiling point, 133
 gas chromatography, 132, 135–136
 infrared spectroscopy, 132
 nuclear magnetic resonance spectroscopy, 132–133
 odor, 132
 polarimetry, 132
 refractive index, 133
 structure of, 129
- Specific unknown, 480
- Spinach
 advance preparation, 145
 carotenoid pigments isolation, 142–148
 chlorophyll isolation of, 142–148
 column chromatography, 144–145
 column running, 146
 fresh vs. frozen, 143
 pigment extraction, 144
 result, analysis of, 147–148
 thin-layer chromatography, 146–147
 TLC plate, 146–147
- Standard-taper ground-glass joints, 618–619
- Steam distillation
 definition of, 802
 immiscible mixtures, 802, 803f
 calculations, 803–804
 methods, 804t
 direct method, 804–806, 805f
 live steam method, 806–807, 806f
 miscible mixtures, 802, 803f
- Steam distillation essential oils, 120–124
 apparatus, 123
 drying, 122
 evaporation, 122
 extraction, 122
 microscale procedure, 121–122
 semimicroscale procedure, 123–124
 yield determination, 122
- Stem correction, 750–752, 752f
- Steric approach control, 291
- Steric energy, 158
- Stone, Edward, 66
- Strain energy, 158
- Sublimation, 797–801
 advantages of, 801
 definition of, 797

methods for, 799–801, 800f
 vacuum sublimation, 799
 vapor-pressure behavior of liquids,
 797–798, 797f
 vapor-pressure behavior of solids, 797–798,
 797f, 798t
 Substitution product, 502
 Substitution reaction, 501
 Sulfa drugs, 423–425
 procedure, 428–429
 sulfanilamide preparation, 426–430
 Sulfur test, 499
 Supernatant liquid, 79
 Surfaces, 174
 Syngas, 249
 Synthetic dyes, 386–391

T

Tautomerism, 177–178
 Tea
 caffeine isolation, 101–103
 analysis, 105–106
 caffeine isolation, 105
 caffeine removal, 105
 solution filtering, 104–105
 solution preparation, 104
 sublimation, 105
 solid phase extraction, 103–106
 Terpenes
 classification of, 117
 history of, 116
 isoprene rule, 117
 isoprene units, 117
 Territorial pheromones, 410
 Theoretical plates, 773
 Thermodynamic product, 457
 Thermoplastics, 432–433
 Thermoset plastics, 433
 Thin-layer chromatography, 88, 828–841
 analgesic drug analysis, 89–93
 commercially prepared TLC plates, 829
 developing TLC plates, 833–834, 834f
 larger thin-layer plates, 831
 microscope slide TLC plates, 830–831
 in organic chemistry, 838–839, 839f
 paper chromatography, 840
 plate development, 47
 plate preparation, 46–47, 47f
 plate spotting, 832f, 833, 833f
 preparative plates, 837
 principles of, 828–829
 reaction, monitoring of, 49–50
 solvent choice, 834–835, 835f
 solvent selection, 48
 thin-layer plate, 828, 829–830
 thin-layer slide, 828, 829–830
 visualization methods, 835–836
 Total energy, 173
 Trail pheromones, 410

Triphenylmethanol, 317–328
 benzoic acid and, 317–328
 procedure, 320–321, 321f, 322–324
 Triphenylpyridine, 366–368
 Michael-Aldol condensation
 reaction, 368
 synthesis of, 368
 Triplet energy transfer, 466

U

Unknown liquids, list of, 65
 Unknowns
 identification of, 480–487
 chemical classification, 481, 482
 chemical classification tests, 484–485
 derivatives, 486
 elemental analysis, 486
 identity confirmation, 486
 inspection, 485
 melting-point determination, 483–484
 optional procedure, 481
 preliminary classification, 483
 preliminary tests, 484
 purification, 484
 solubility behavior, 484
 spectroscopy, 481, 482, 485–486
 solubility tests, 487–493
 Unleaded gasoline, 237
 Unsaturation tests, 500–503
 simple multiple bonds, 501–503

V

Vacuum distillation, 785–796
 bulb-to-bulb distillation, 791–792, 791f
 closed-end manometer, 793–794, 793f, 794f
 macroscale equipment, 788–790, 789f
 manometer use, 795–796, 795f
 mechanical vacuum pump, 792–793, 792f
 microscale methods, 785–788, 786f, 787f
 microscale vacuum distillation, 788–791
 rotary fraction collectors, 791, 791f
 semimicroscale equipment, 788–790
 simple microscale apparatus, 788
 Vane, J.R., 67
 Vanillin
 esterification reactions of, 586–588
 procedure, 587
 Variation Principle, 168
 Vermifuge, 73
 Vision, chemistry of, 137–140
 Visualization method, 829, 835–836
 VOCs. *See* Volatile organic compounds
 Volatile organic compounds (VOCs), 239
 Volume measurement, 632–639
 automatic pipettes, 633, 633f
 balances, 638–639, 638f
 beakers, 638
 conical vials, 638
 dispensing pumps, 633–634, 634f

Volume measurement (*Continued*)

- graduated cylinders, 638
- graduated pipettes, 634–637, 635f, 636f
- Pasteur pipettes, 636–637, 637f

Volume measurement (*Continued*)

- pipette pump, 635f
- syringes, 637

Vomeronasal organ, 411

Vooshall, Leslie, 414

W

Wald, George, 137

Water bath, 3

Water-jacketed condenser, 649

Wavefunction, 169

Wavenumbers, 880

Weight measurements, 632–639

Wittig reaction, 369

Wittig salt, 372

Wort, 149–150

X

Xanthophylls, 142

Y

Ylide, 306, 370

Z

Zeiss polarimeter, 870, 871f

Zymase, 152

This is an electronic version of the print textbook. Due to electronic rights restrictions, some third party content may be suppressed. Editorial review has deemed that any suppressed content does not materially affect the overall learning experience. The publisher reserves the right to remove content from this title at any time if subsequent rights restrictions require it. For valuable information on pricing, previous editions, changes to current editions, and alternate formats, please visit www.cengage.com/highered to search by ISBN, author, title, or keyword for materials in your areas of interest.

Important notice: Media content referenced within the product description or the product text may not be available in the eBook version.